

# Microdeflectometry – a novel tool to acquire 3D microtopology with nanometer height resolution

Gerd Häusler,\* Claus Richter, Karl-Heinz Leitz, and Markus C. Knauer

*Institute of Optics, Information and Photonics, Max Planck Research Group,  
University Erlangen-Nuremberg, Staudtstr. 7/B2, 91058 Erlangen, Germany*

*\*Corresponding author: haeusler@physik.uni-erlangen.de*

Compiled March 31, 2022

We introduce "microdeflectometry", a novel technique for measuring the microtopology of specular surfaces. The primary data is the local slope of the surface under test. Measuring the slope instead of the height implies high information efficiency and extreme sensitivity to local shape irregularities. The lateral resolution can be better than one micron whereas the resulting height resolution is in the range of one nanometer. Microdeflectometry can be supplemented by methods to expand the depth of field, with the potential to provide quantitative 3D imaging with SEM-like features.

*OCIS codes:* 120.0120, 120.3940, 120.6650, 180.0180, 180.6900.

We introduce a microscopic adaptation of deflectometry [1]. Specifically, we modify "Phase-Measuring Deflectometry" (PMD)[2]. "Microdeflectometry" provides quantitative slope images with lateral resolution better than one micrometer and slope resolution in the range of 1 mrad. A surface height variation of one nanometer can be detected within the diffraction-limited resolution cell of the optical system. The method is incoherent (low noise) and provides an angular dynamical range as big as the aperture of the microscope.

PMD and similar methods [3, 4, 5] is based on the observation of mirror images of remote patterns, using the object under test as a mirror. The mirrored patterns are distorted, depending on the shape of the object. If we use sinusoidal fringes and a calibrated system, the local slope of the surface can be calculated at about one million object points (x,y), within a few seconds. To get the best possible lateral resolution we focus on the surface of the object under test. Due to the limited depth of field (dof) we cannot focus on the object and on the screen at the same time. However, the measured phase will not be altered by the blurring since the fringes are sinusoidal [2, 6].

Measuring the local slope instead of the height of an object has significant advantages. Optical range sensors commonly need a great deal of 2D raw data to finally calculate the 3D image [7]. They are inefficient in terms of information theory. In contrast, deflectometric systems perform significant source encoding and display high information efficiency. Moreover, local slope data are more effective in detecting local details and defects than height data [8]. Macroscopic PMD is a now-established technique for measuring macroscopically specular free-form surfaces, specifically eye glass lenses [9]. Other applications are painted car bodies, car window glasses or precision-machined surfaces.

Our aim is to adapt this measuring principle to a microscopic field of view down to 100  $\mu\text{m}$ . This opens up new applications. A few examples are shown in

this paper: Cutting tools (Fig. 3), microoptical elements (Fig. 4), and wafers (Fig. 5). With high aperture and high resolution, many surfaces that are matt in a macroscopic setup display specular reflection, and these can be measured with microdeflectometry. With high-aperture imaging steep slopes are accessible (up to  $\pm 60^\circ$ ).

We now explain how to implement microdeflectometry. It is not possible to use a diffusing remote screen, as in macroscopic PMD, to display the fringe pattern, because of the small working distance of a microscope. Our solution is different: As depicted in Fig. 1, we project an aerial image at a distance  $d$  from the object. This projection is achieved by means of the same microobjective that acquires both the image of the object and the mirror image of the fringe pattern. We generate this pattern using an electronically controllable light modulator and a beam splitter, neither shown in Fig. 1. This scheme enables an incoherent (low optical noise) and coaxial (no shadowing) illumination with high aperture (high resolution, big angular range) and small working distance.

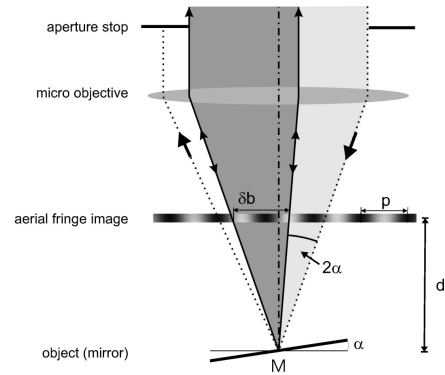


Fig. 1. Schematic setup for microscopic PMD

We briefly explain the imaging scheme for the example of a point  $M$  located at the object surface. The cone of

imaging rays is represented by the dark gray area. The observed image intensity is given by the intensity within this cone, averaged over the distance  $\delta b$ . The tilted object causes vignetting and the rays within the light gray cone do not contribute to imaging. The cone of imaging rays depends on the slope  $\alpha$ . Hence the observed phase  $\varphi$  depends on the slope, with  $\varphi = 2\pi d\alpha/p$ . Note that due to the vignetting, microdeflectometry is only half as sensitive as macroscopic deflectometry.

We now discuss the physical limitations. In [2, 6] an uncertainty relation as given in Eq. (1) was derived that connects the lateral resolution  $\delta x$  and the angular measuring uncertainty  $\delta\alpha$ , where  $\lambda$  is the wavelength and  $Q$  is a quality factor, depending solely on the camera noise:

$$\delta x \cdot \delta\alpha = \lambda/Q \quad (1)$$

Equation (1) appears a little disappointing since we learn that with higher lateral resolution (with higher aperture) the slope uncertainty will become bigger. However, if we put in the experimental value of  $Q \approx 1000$ , we get approximately:

$$\delta x \cdot \delta\alpha = 1\text{nm} \quad (2)$$

Equation (2) reveals a remarkable result: The minimum height variation  $\Delta z$  that can be detected within the diffraction-limited resolution distance of the microscope is constant and in the range of only 1 nm. One might doubt that this relation is still valid for a high-aperture microscope, since the dof may be below  $1\ \mu\text{m}$ , and the distance  $d$  might become bigger than a hundred times the dof. Yet the result of our experiments, shown in Fig. 2, demonstrates that Equation (2) is valid even for very high aperture and a very small dof. This is possible because we scale up the fringe period  $p$  properly with the distance  $d$  in order to get sufficient contrast in the fringe image. (It should be mentioned that there are systems for pointwise macroscopic measurement [11] with  $Q > 10^5$ .)

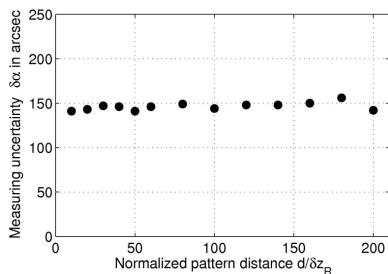


Fig. 2. Experimental measuring uncertainty  $\delta\alpha$  vs. normalized pattern distance  $d$ .  $\delta z_R$  is the Rayleigh depth of field. Objective: 50x0.85

Now we present some examples of intensity-encoded slope images. The first example is a razor blade (Fig. 3). The measuring field is 1.4 mm (objective 5x0.12). Figure 3 a displays the different slope angles of the blade, in y

direction. The measured angles are close to the reference values of the manufacturer ( $2.0^\circ$ ,  $4.5^\circ$  and  $6.0^\circ$ ). In Fig. 3 b the slope in x direction is shown, with machining traces emphasized.

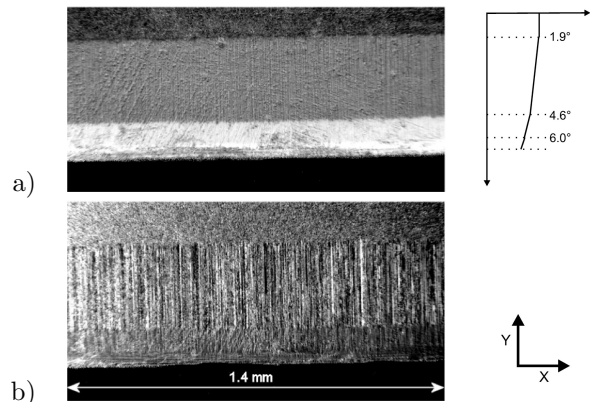


Fig. 3. Razor blade. a) Local slope in y direction, b) in x direction.

The second example is an array of cylindrical microlenses (Fig. 4 a). The measuring field is  $105\ \mu\text{m}$  (objective 50x0.85). The edges of the lenses have a slope of about  $11^\circ$ . The scratches (see arrow) are about 40 nm deep. In Fig. 4 b a cross section of the height map obtained by (two-dimensional) integration of the slope data [10] is depicted. The profile data are degraded (error  $< 10\%$ ) by the still-preliminary calibration of the system. On the right side of the profile there is a systematic error due to defocusing. We later will show that the defocusing problem can be overcome.

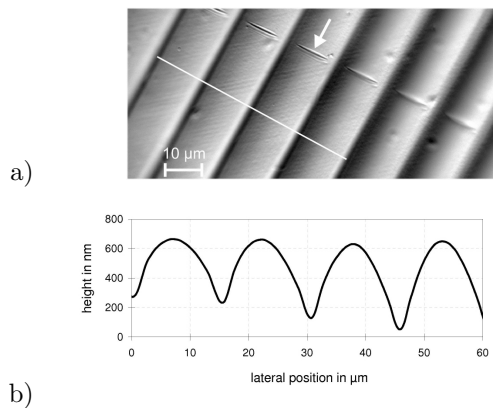


Fig. 4. a) Slope image of a microlens array; b) cross section through height map (by integration).

As a further example, Fig. 5 displays the topology of a wafer.

Further on this paper, we will briefly address the problem of the limited depth of field. Figure 6 demonstrates that it is possible to combine microdeflectometry with the expansion of dof. The result is an SEM-like image,

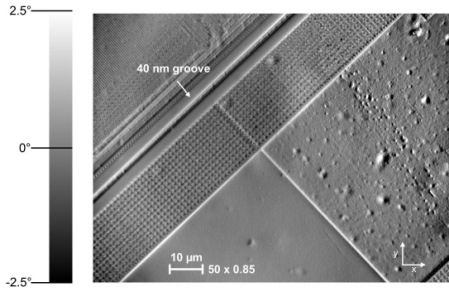


Fig. 5. Wafer section. Local slope in y direction.

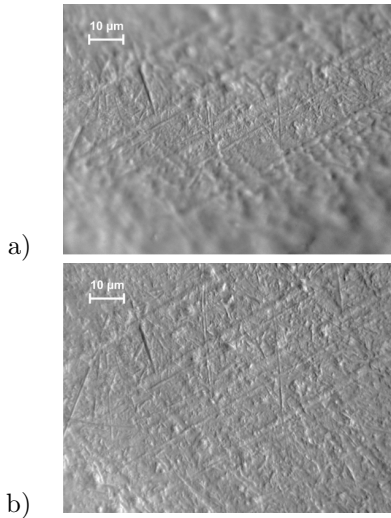


Fig. 6. Microdeflectometry with additional expansion of depth of field. Slope image (x direction) of a used bearing ball. Height range approx.  $17\ \mu\text{m}$  (30-fold dof, objective  $50\times 0.85$ ). To reduce the dynamical range of the picture, only the deviation of the measured slope from the slope of the sphere is displayed. a) standard measurement; b) Measurement with about 30-fold extended dof.

which additionally provides quantitative 3D data. We will discuss this method in detail in an additional paper.

We conclude that microdeflectometry displays several interesting features:

a) The inherently measured signal is the slope. Slope encoding makes DIC images (differential interference contrast [12]) and SEM images (scanning electron microscope) so intriguing. The visual impression is significantly advantageous over a pure intensity image. The slope encoding enhances extremely small depth variations of the surface, down to the nanometer and even sub-nanometer regime. From the slope, the shape  $z(x,y)$  can be calculated with extremely low noise [10], since integration is a noise-suppressing operation. (We admit that steps at an object surface cannot be measured. Most smooth surfaces do not display steps, fortunately). High-accuracy quantitative microdeflectometry still requires an improved calibration procedure.

b) The noise level is small. Microdeflectometry can be implemented with high illumination aperture, so resid-

ual coherent noise artifacts will have very low contrast. Moreover, the high aperture enables a high dynamical range of the slope measurement. In our experiments we achieved a dynamical range of up to 2500. This is much better than what can be achieved in intensity images.

c) The lateral resolution and the angular measuring uncertainty are coupled via an uncertainty relation, so that the corresponding parameters can be adapted as required by the application. For eye glass measurements the accuracy of the refractive power ( $\sim$ curvature) is more valuable than the lateral resolution. For semiconductor applications, the lateral resolution is more important.

d) Depth of field: the small dof is a major drawback of microscopic imaging. It is intriguing and a challenge to create SEM-like images (big dof) *optically* with microdeflectometry. We have demonstrated that this is possible, by combining microdeflectometry with the expansion of the dof. In fact, the expansion of the dof with microdeflectometry is easier than with conventional microscopy.

Microdeflectometry provides appealing pictures, with high lateral resolution, nanometer sensitivity for local surface features, low noise and quantitative 3D features. The technology is simple, incoherent, and has the potential to compete against interferometry.

## References

1. R. Ritter, R. Hahn, "Contribution to analysis of the reflection grating method", in *Optics and Lasers in Engineering*, **No.1**, pp 33-44 (1983).
2. G. Häusler, German patent DE 19944354.
3. D. Pérard, J. Beyerer, "Three-dimensional measurement of specular free-form surfaces with a structured-lighting reflection technique", *Proc. SPIE* **3204**, 74-80 (1997).
4. M. Petz, R. Tutsch, "Reflecting grating photogrammetry", in *VDI-Berichte No. 1844*, pp. 327-338 (2004).
5. T. Bothe, W. Li, C. von Kopylow, W. Jüptner, "High Resolution 3D Shape Measurement on Specular Surfaces by Fringe Reflection", *Proc. SPIE* **5457**, 411-422 (2004).
6. M. C. Knauer, J. Kaminski, G. Häusler, "Phase Measuring Deflectometry: a new approach to measure specular free-form surfaces", *Proc. SPIE* **5457**, 366-376 (2004).
7. C. Wagner, G. Häusler, "Information theoretical optimization for optical range sensors", *Appl. Opt.*, **42 No. 27**, 5418 - 5426 (2003).
8. M. C. Knauer, T. Bothe, S. Lowitzsch, W. Jüptner, G. Häusler, "Höhe Neigung oder Krümmung?", in *DGaO-Proc.*, **107**, B30 (2006).
9. 3D-Shape GmbH, [http://www.3d-shape.com/produkte/pmd\\_e.php](http://www.3d-shape.com/produkte/pmd_e.php)
10. S. Lowitzsch, J. Kaminski, M. C. Knauer, G. Häusler, "Vision and Modeling of Specular Surfaces", *Proc. VMV*: 479-486 (2005).
11. R. D. Geckeler, I. Weingärtner, "Sub-nm topography measurement by deflectometry: flatness standard and wafer nanotopography", *Proc. SPIE* **4779**, 1-12 (2002).
12. M. Davidson, and M. Abramowitz, "Optical microscopy" in *Encyclopedia of Imaging Science and Technology*, J. Hornak (ed), Wiley, N. Y., 1106-1141 (2002).