

Companding technique for high dynamic range measurements using gafchromic films.

Frank Van den Heuvel,¹ Wouter Crijns,¹ and Gilles Defraene¹

*Department of Experimental Radiotherapy,
Catholic University of Leuven,
Leuven, Belgium*

We propose a methodology to perform dose measurements using gafchromic films which can span several decades of dose levels. The technique is based on a rescaling approach using different films irradiated at different dose levels. This is combined with a registration protocol correcting positioning and scaling factors for each film.

The methodology is validated using TLDs for out-of-field doses. Furthermore, an example is provided using the technique to characterize a small 7.5 mm sized radiosurgery cone and compared to measurements made with a pinpoint chamber.

I. INTRODUCTION

The use of gafchromic films as a measurement device has become popular during the last few years. There are many reasons for this among which: 1) high resolution, 2) no additional development needed, 3) off the shelf reading tools, and 4) tissue equivalence. Some disadvantages include: 1) variations in sensitivity 2) heterogeneity within a film does exist, and 3) a limited dynamic range does not allow the films to be used to quantify beam characteristics or to measure relatively low doses. Of these disadvantages the last is the most important as batch production and pre-scanning allow to minimize the first two problems.

In this paper we describe a technique to increase the dynamic range to low doses by making use of a rescaling algorithm inspired by companding noise reduction. Companded transmission is a technique that has been used to transmit data over noisy telegraph lines³ and has also been used in the audio-industry to reduce cassette tape noise⁶. In essence a signal or part of a signal is amplified before transmission after which it is reduced to its normal level before being fed to the output circuit.

In addition, we will show that this approach allows to perform small field measurements with superior resolution than heretofore possible.

II. METHODS AND MATERIALS

A. Films

We use EBT-II gafchromic films provided by International Specialty Products¹. Scanning the films at a resolution of 150dpi using a 10000XL Epson scanner. The films are batch calibrated to provide absolute dose values. Two films represent the complete batch and is calibrated using a two film methodology. This procedure has been described elsewhere⁵ and provided an absolute error level of the order of 3 to 5 using the red channel output%. The films are calibrated to provide accurate dosimetry between 0.2Gy and 4 Gy.

B. Companding technique

In order to provide correct data over several decades of dose (i.e. from 4 to 0.002Gy) we use different films to and expose them to different levels of dose. The reference film's dose

(D_1) is chosen in such a way that the maximal expected dose is close to the maximum from the dose range $[U_1, L_1]$ for which the film was calibrated, with U_1 being the upper limit and L_1 the lower limit. The subsequent films are irradiated with doses that are higher than the previous film using a multiplication factor F_i . Where we denote the dose for the i -th film: $D_i = F_{i-1}D_{i-1}$. The choice of this factor is made in such a way that there is information overlap between the film. The emphasis factor (F_i) used for the $i + 1$ -th film should be chosen so that $L_i < U_i/F_i$.

All films are converted to dose using a calibration protocol that links the transmittance to dose. The specifics of this calibration procedure are extensively dealt with elsewhere⁵.

C. Registration

As the films exhibit a high resolution (1 pixel \sim 0.16mm) it is difficult to exactly reproduce the position of a given film. Additionally, variations in the beam output of a linear accelerator introduces an uncertainty in the dose delivered and thus in the factor used to de-emphasize the signal. Finally, variations in film sensitivity within a single batch of films have been shown to exist and need to be taken into account. For the previously mentioned reasons it might occur that perfect alignment and dose levels differ from film to film. To solve this we introduced an alignment protocol taking into account in both dosimetric and spatial variations. This is not a straightforward problem as the geometric registration and the dose level are not independent parameters.

The registration between two films of subsequent level of companding is performed as follows:

- We define the registration as a 4-dimensional vector $(X_i, Y_i, \theta_i, F_i)$ with X_i and Y_i being the translational coordinates and θ_i a rotational parameter, the component F_i denotes the emphasis factor.
- The films are converted to dose using the aforementioned calibration. Only data that falls within the upper limits (U_i) and lower limit (L_i) of calibration is retained.
- The two films are registered using data that is contained in the overlap between the two calibration ranges. This is dependent on the position of the film as well as the dose factor.

- The cost function to register the images is the dose difference between the two images multiplied by the gradient image $G = 1 + a \times \sqrt{G_X^2 + G_Y^2}$, with G_X and G_Y the gradient calculated in the X and Y direction of the film which are chosen to follow the film edges. The parameter a is set at 500 to allow the gradient function to have an impact on the images used as the gradient image is normalized to have an integrated amplitude of 1.
- The cost function is optimized using a constrained minimization algorithm as implemented in the function *fmincon* in Matlab⁴, with the option *interior-point approach*.
- The minimization procedure is performed in a progressive manner starting from sub-sampled image at 50dpi, over 75dpi and finally a full resolution image at 150dpi. For every increase in resolution the constraints within which the optimization is performed are tightened. They are respectively (5, 5, 10°, 10%), (0.1, 0.1, 0.5°, 0.5%), and (0.03, 0.03, 0.05°, 0.15%). Using the same notation as for the transformation vector with all spatial values in cm.

The end result is a dose image of high resolution and dynamic range that can be used to characterize a radiation beams in treatment planning systems or to provide data on radiation dose levels far from the field edge. An example of a resulting data set is provided in Figure 1.

D. Validation

The validation is performed with a single set-up. A stack of solid water is irradiated using a clinical 6MV photon beam from a Varian C/D 2100 linear accelerator. The gafchromic films are placed at a depth of 5cm, with the surface of the stack located at 95cm. We present three different tests.

1. A 10×10 cm² field is compared to the results from a planning system (Varian EclipseTM) with the pencil beam and the AAA-algorithm. In addition the low dose area is compared to measurements using TLDs. This data set has also been used in the Allegro project to compare with TPS data. Here we add the TLD validation.

2. A BrainlabTM radio surgery cone with an opening size of 30mm with the machine jaws set to a $5 \times 5 \text{ cm}^2$ size is compared to the data obtained for entry in the BrainscanTM planning system. For this purpose the field is measured with a PTW pinpoint chamber type 31006 of size 0.015 cc^2 using a waterphantom.
3. In an identical setup a cone of opening size 7.5mm was also assessed.

1. TLD

LiF:Mg,Ti TLDs (TLD-100 EXTRAD, Harshaw) are used. Calibration of the dosimeters was performed in a 6 MV beam at 5cm depth in a polystyrene phantom. The TLDs are read-out in a Harshaw 6600 hot gas reader, equipped with a $^{90}\text{Y}/^{90}\text{Sr}$ source. In order to minimize initial fading effects, read-out of calibration TLDs that have been irradiated with the internal source simultaneously with the experiment, is part of the standard processing procedure.

III. RESULTS

A. 10×10 fields

Figure 2 provides a graph for the open field test. It is clear that the TLD's agree well with the film measurements. Due to the low dose on the TLD's there are considerable absolute error bars on these measurements. Also with doses to TLDs of this magnitude the contribution of the background can cause an offset in the measurements. A good agreement with the planning system is found when AAA is used. The pencil beam data does not agree that well.

B. Small field dosimetry 30mm cone

As in the classical field one can notice the difference in resolution of the gafchromic film compared to the measurement. In this case a pinpoint chamber was used. Figure 3 shows a log plot of the relative dose. The measurement here is performed using a pinpoint chamber, however a broadening of the field edge compared to the film is easily discernible.

C. Small field dosimetry 7.5mm cone

In the case of the data obtained with irradiation using the smallest cone there is no plateau in the data measured using the pinpoint chamber, leading us to believe that there are partial volume effects in play already. The film measurement does show a plateau phase. Also the ratio of the maximal measured dose compared to the output from a 10x10 field (output ratio) was more than double with film than with the pinpoint chamber. Due to these differences we have never used the smallest cone clinically. Furthermore, away from the central axis we note an increase in dose. Although the dose remains low increasing from 0.2% to 0.5% of the dose delivered to the central axis.

D. Alignment

In table I the transformation needed to match the different films show that there was a difference in response or machine output of the films used to measure the smaller dose compared to the normal film. It also shows that reasonable accuracy can be reached by mechanically positioning the films correctly in treatment as well as readout position. However, although shifts are of the order of less than 0.5 mm in the case of very steep gradients important inconsistencies can be introduced when using manual alignment only. The latter becomes important in the case of very small fields.

IV. DISCUSSION

The approach as outlined here has been used to perform measurements within the framework of the ALLEGRO project. The goal of one of the workgroups was to determine the dose outside of the high dose volume. This to estimate the impact of newer radiotherapy techniques on the whole body dose. A second question to answer was whether data from commercial planning systems would be able to serve as a resource to obtain information on low-dose exposures. It was therefore necessary to enhance the accuracy of low dose measurements.

From the work in this paper and the comparisons with TLD measurements and ion-chamber data, we conclude that this technique provides data that not only allows estimates of low-dose contributions to patients, but also provides spatial information. The use of

gafchromic film in this fashion will also allow its use in the characterization of radiation systems that deliver doses to very small volumes.

ACKNOWLEDGEMENTS

This study was partly funded by the ALLEGRO-project (7th Euratom Framework Programme for Nuclear Research and Training, Fission-2008-3.2.1).

AUTHOR CONTRIBUTIONS

Mr. Van den Heuvel outlined the idea and approach and provided a proof of concept. Mr Crijs actualized the registration component. Mr Defraene provided the TLD validation and the registration test. All authors contributed to discussions and writing the article.

REFERENCES

¹ISP, New Jersey, USA.

²PTW Freiburg, Germany.

³Clark, A. (1928). Electrical picture-transmitting system. *United States patent*, 1691147.

⁴Coleman, T. F. and Li, Y. (1994). On the convergence of interior-reflective Newton methods for nonlinear minimization subject to bounds. *Mathematical Programming*, 67(1):189–224.

⁵Crijs, W. and Van den Heuvel, F. (2011). Gafchromic film dosimetry: calibration methodology and error analysis. *Submitted to Med Phys*.

⁶Dolby, R. (1967). An audio noise reduction system. *J. Audio Eng. Soc*, 15(4):383–388.

TABLES

10×10 Field	Decade 1–2	Decade 1–3
X–shift	0.65 mm	1.99 mm
Y–shift	0.64 mm	3.35 mm
Rotation	-0.89°	-1.29°
Correction for F_i	1.037	1.071
Cone 30mm	Decade 1–2	Decade 1–3
X–shift	-0.00 mm	-0.24 mm
Y–shift	-0.15 mm	0.36 mm
Rotation	-0.04°	-0.33°
Correction for F_i	0.934	0.957
Cone 7.5mm	Decade 1–2	Decade 1–3
X–shift	0.58 mm	-0.45 mm
Y–shift	0.07 mm	1.53 mm
Rotation	-2.99°	-2.76°
Correction for F_i	1.08	1.17

TABLE I. Transformations performed to align the different films in this setup. The third film is matched to the second film. The data given are for all films with respect to the first film. This implies that errors for the third film are cumulative.

FIGURE CAPTIONS

FIG. 1. Combination of three decades of irradiation data for a 10×10 cm² sized beam. The z-axis is logarithmic and in cGy, the other axes are in cm

FIG. 2. Validation of a classic radiation therapy beam with a 6MV energy. The AAA implementation follows the measurement excellently. Moreover the TLD measurements show good agreement. A slight offset can be expected due to some background dose built up in the TLDs. Also note the difference in the field edge between planning and film measurement. The TPS data was modelled using data obtained with a compact IC 10 ion chamber (Wellhöfer).

FIG. 3. A 30mm cone delivers 6MV photons on the films. We compare to the results from a pinpoint chamber. The blue, orange and black squares delineate the spatial and dosimetric limitations for respectively the first, second, and third decade. Note that even while the pin-point chamber is small a difference can be seen in the resolving power compared to film measurements.

FIG. 4. A 7.5mm diameter cone delivers 6MV photons on the films. We compare to the results from a pinpoint chamber. The blue, orange and black squares delineate the spatial and dosimetric limitations for respectively the first, second, and third decade. Note that even while the pin-point chamber is small a difference can be seen in the resolving power compared to film measurements. Also note that further away from the central axis in the region not measured by the ion chamber the dose increases by an order of magnitude.

FIGURES

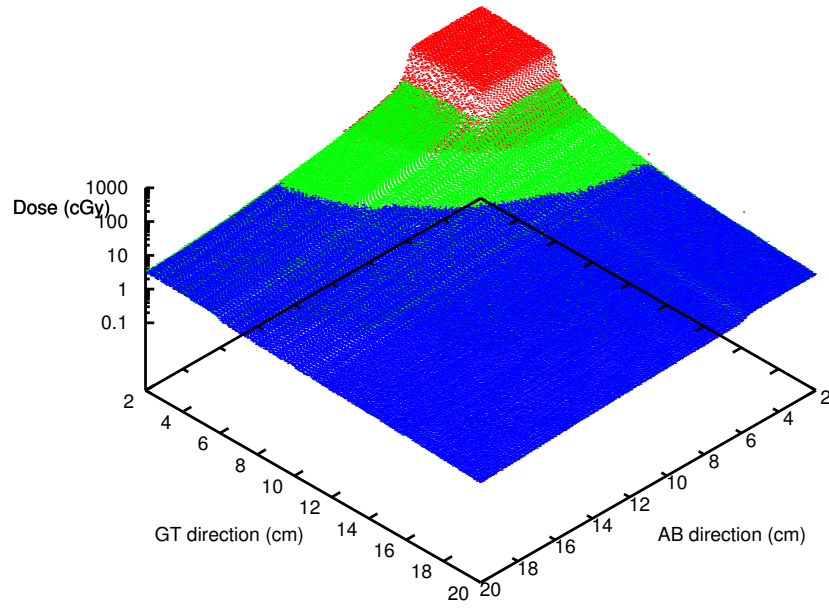


FIG. 1. Combination of three decades of irradiation data for a $10 \times 10 \text{ cm}^2$ sized beam. The z-axis is logarithmic and in cGy, the other axes are in cm

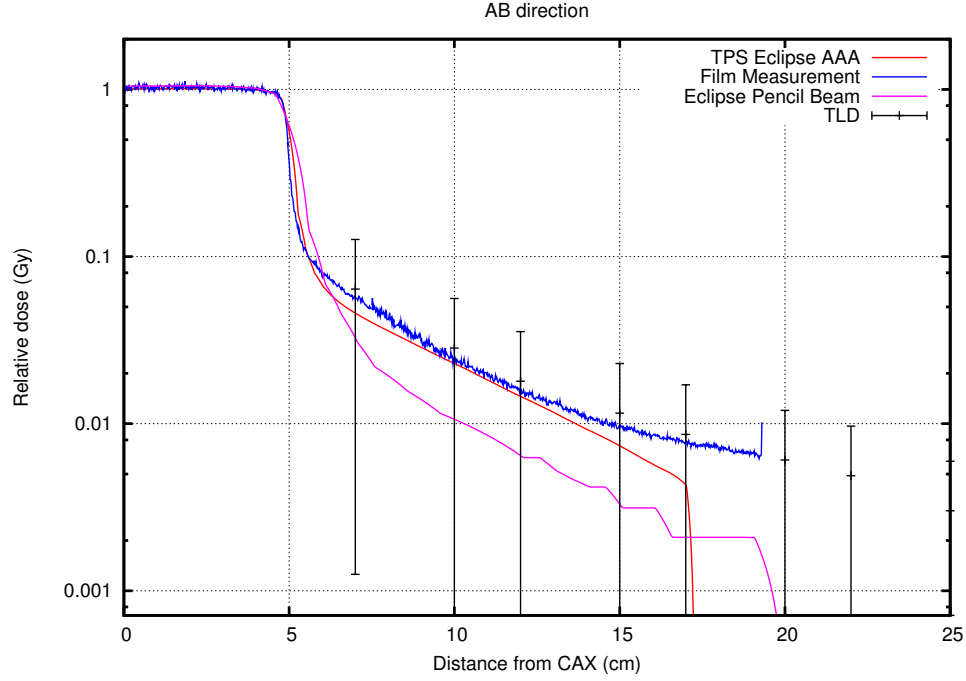


FIG. 2. Validation of a classic radiation therapy beam with a 6MV energy. The AAA implementation follows the measurement excellently. Moreover the TLD measurements show good agreement. A slight offset can be expected due to some background dose built up in the TLDs. Also note the difference in the field edge between planning and film measurement. The TPS data was modelled using data obtained with a compact IC 10 ion chamber (Wellhöfer).

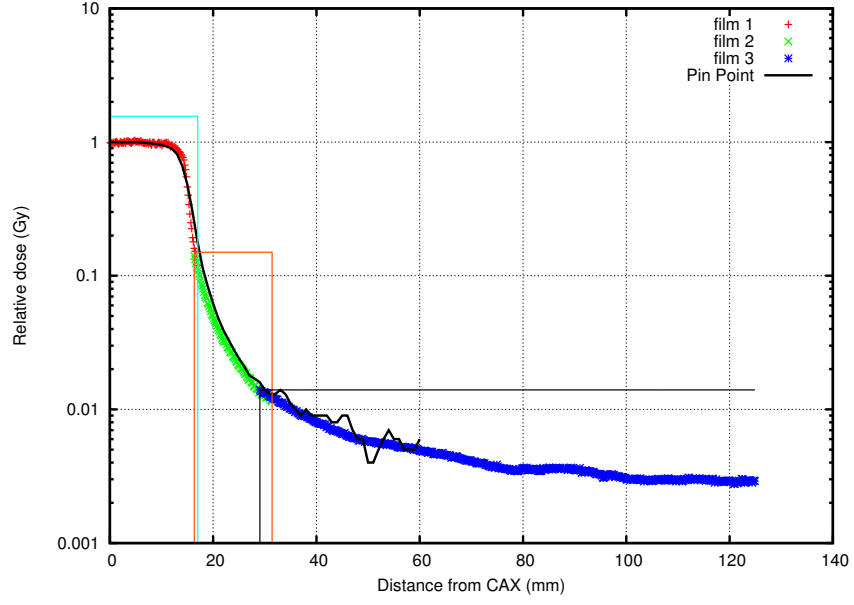


FIG. 3. A 30mm cone delivers 6MV photons on the films. We compare to the results from a pinpoint chamber. The blue, orange and black squares delineate the spatial and dosimetric limitations for respectively the first, second, and third decade. Note that even while the pin-point chamber is small a difference can be seen in the resolving power compared to film measurements.

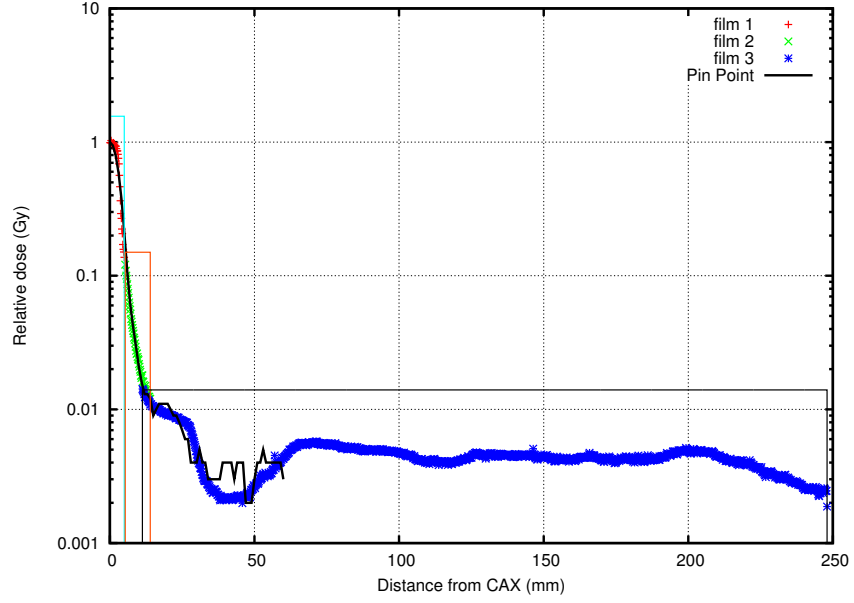


FIG. 4. A 7.5mm diameter cone delivers 6MV photons on the films. We compare to the results from a pinpoint chamber. The blue, orange and black squares delineate the spatial and dosimetric limitations for respectively the first, second, and third decade. Note that even while the pin-point chamber is small a difference can be seen in the resolving power compared to film measurements. Also note that further away from the central axis in the region not measured by the ion chamber the dose increases by an order of magnitude.