1 Cluster SIMS Microscope Mode Mass Spectrometry 2 Imaging

3 András Kiss¹, Donald F. Smith¹, Julia H. Jungmann¹, Ron M.A. Heeren^{1*}

4 ¹ FOM Institute AMOLF, Science Park 104, 1098 XG Amsterdam, The 5 Netherlands

6 * Author to whom correspondence should be addressed. Email: heeren@amolf.nl

7

8 Abstract

9 Rationale

Microscope mode imaging for secondary ion mass spectrometry is 10 a technique with the promise of simultaneous high spatial resolution and 11 12 high speed imaging of biomolecules from complex surfaces. Technological developments such as new position-sensitive detectors, in combination 13 with polyatomic primary ion sources, are required to exploit the full 14 potential of microscope mode mass spectrometry imaging, i.e. to efficiently 15 push the limits of ultra-high spatial resolution, sample throughput and 16 17 sensitivity.

18 Methods

In this work, a C₆₀ primary source is combined with a commercial mass microscope for microscope mode secondary ion mass spectrometry imaging. The detector setup is a pixelated detector from the Medipix/Timepix family with high-voltage post-acceleration capabilities. The system's mass spectral and imaging performance is tested with various benchmark samples and thin tissue sections.

25 **Results**

We show that the high secondary ion yield (with respect to "traditional" monatomic primary ion sources) of the C_{60} primary ion source and the increased sensitivity of the high voltage detector setup improve microscope mode secondary ion mass spectrometry imaging. The analysis time and the signal-to-noise ratio are improved compared to other microscope mode imaging systems, all at high spatial resolution.

32 Conclusions

We have demonstrated the unique capabilities of a C₆₀ ion microscope with a Timepix detector for high spatial resolution microscope mode secondary ion mass spectrometry imaging.

1 Introduction

Mass spectrometry imaging (MSI)^[1, 2] is a technique for the label 2 free study and visualization of the distribution of multiple molecular 3 species on complex surfaces, such as thin tissue sections. It has two main 4 advantages over other common imaging techniques used. First it has 5 chemical identification capabilities and second, no a priori knowledge of 6 the sample is required. Matrix assisted laser desorption ionization (MALDI) 7 has seen widespread use for mapping of intact biological molecules from 8 complex surfaces. Spatial resolution in such experiments is typically 9 limited to 10-50 µm, due to laser beam focusing. However, secondary ion 10 mass spectrometry (SIMS) has a long history of sub-micrometer chemical 11 imaging of a variety of sample substrates^[3-9]. 12

SIMS is the oldest ionization method used for mass spectrometry 13 imaging^[10]. Early primary ion sources for SIMS used atomic primary ions 14 such as Ga⁺, Cs⁺ or In⁺. These ion sources have the advantage that the 15 ion beam can be focused to a very small spot size, thus a very high spatial 16 resolution can be achieved (10s of nanometers)^[11, 12] for MS imaging 17 experiments. A drawback of these sources is the high degree of 18 fragmentation of the secondary ions. The range of detectable ions is 19 limited to elemental ions and small organic fragments such as CH_3^+ or CN^- 20 , albeit at high spatial resolution. For many years, this restricted the use of 21 22 SIMS to solid state physics and to the study of semiconductor surfaces. The introduction of polyatomic primary ion sources was one of the biggest 23 advancements in the field^[13-15]. It was shown by several groups that these 24 sources have a higher secondary ion yield and they provide "softer" 25 ionization^[13, 16-20], thus opening the field of SIMS for biological research 26 applications. However, these ion sources are more difficult to focus due to 27 the strong space-charge effect associated with the combination of high 28 primary ion current and small spot size. Thus, a compromise is needed 29 between the primary ion current which is related to the secondary ion yield 30 and spatial resolution. This is shown by the fact that the highest spatial 31 resolution achieved with a C_{60} primary ion source was 1 µm as opposed to 32 the tens of nanometers with monoatomic primary ion beams. This was 33 demonstrated by the Vickerman group which built a C₆₀ buncher Time-of-34 Flight (TOF) instrument for high spatial resolution imaging^[19]. 35

Typically, mass spectrometry imaging is performed in the microprobe mode, where a focused laser (MALDI) or primary ion beam (SIMS) measures the sample pixel by pixel. The spatial resolution is dependent on the spot size of the ion source in microprobe mode. *Microscope* mode imaging does not require highly focused beams. Rather, large beam sizes (usually around 200-300 μ m) are used to desorb and ionize the molecules on the sample surface. After the desorption and

ionization event, the spatial distribution of the ions is preserved as they 1 travel through the mass spectrometer and then the ion image is projected 2 onto a position-sensitive detector. The main advantage of this approach is 3 that the spatial resolution is independent from the size of the laser focus or 4 primary ion beam. In this way the analysis speed is significantly increased 5 6 for the same area, spatial resolution and repetition rate compared to 7 microprobe mode imaging. Additionally, it eases the difficulty of focusing the laser beam or primary ion source to very small spot sizes (few µm or 8 less). 9

So far, the widespread use of microscope mode mass spectrometry 10 imaging was limited by the lack of an appropriate position-sensitive 11 12 detector with simultaneous measurement of time-of-flight and position information. The first detector used for microscope mode imaging 13 consisted of the combination of a dual microchannel plate (MCP) with a 14 phosphor screen and a charge-couple-device (CCD) camera^[21, 22]. The 15 main limitation of these detector systems is that they are not capable of 16 recording the time-of-flight of the arriving ions, only their position. Thus, 17 the selection of an ion of interest with an electrostatic blanker in the mass 18 spectrometer is needed to record a selected-ion image. Because of the 19 sample damage due to the ionization process, this is only possible for a 20 few ions before the sample is depleted. It is also very time consuming 21 22 because the sample needs to be imaged separately for every ion of interest. 23

Recently, other detectors, such as delay-line detectors, have been 24 tested for microscope mode imaging^[23, 24]. These detectors have the 25 capability to record the time-of-flight and the spatial position of an ion 26 27 simultaneously. However, they lack multi-hit capabilities for typical mass spectrometry imaging event rates. As a result, these detection systems 28 require low ion loads and are only well suited for SIMS experiments and 29 not MALDI. Additionally, there is no direct feedback during measurement 30 because of the time consuming image reconstruction process which 31 makes optimization of measurement parameters difficult. 32

The latest development in the field of microscope mode mass 33 spectrometry imaging has been the introduction of pixelated detectors, 34 such as the Medipix/Timepix detector family^[25-30]. In these detectors, every 35 pixel acts as an individual detector capable of recording the time-of-flight 36 of the arriving ions with respect to an external trigger signal. Combined 37 with MCPs, they offer multiplexed ion detection capabilities where every 38 ion hit is registered by multiple pixels. This multiplexed detection results in 39 increased sensitivity for these systems. Additionally, the spatial 40 information is determined by the pixel address of each pixel. The 41

capabilities of the Timepix detectors were previously demonstrated for
 microscope mode MALDI imaging^[27-30].

In earlier work, we presented the first example of microscope mode 3 SIMS imaging with a Timepix detector^[31]. The potential of this detection 4 system for SIMS microscope mode MSI, namely superior spatial resolving 5 power and signal-to-noise was demonstrated. Also, these first experiments 6 revealed that this initial setup lacked the sensitivity and speed needed for 7 practical use in biological studies. Limitations were the low secondary ion 8 yield of the gold primary ion source, the lack of ion post-acceleration in the 9 Timepix setup and the low repetition rate of the Timepix detector. Also, 10 negative mode measurements were not possible with this setup. 11

This work presents an improved Timepix detector setup^[29] with ion 12 post-acceleration and negative ion detection capabilities in combination 13 with a C₆₀ primary ion source. This combination offers improvements in 14 most of the crucial parts of the system such as sensitivity, speed and 15 mass range. In particular, the use of a polyatomic C_{60} primary ion source 16 promises improved ion yield performance as compared to the gold liquid 17 18 metal ion source and is evaluated in detail in this work. The capabilities of the system for microscope mode C₆₀ SIMS imaging were demonstrated 19 with various benchmarks and thin tissue sections. 20

21 Materials and methods

22 The Medipix/Timepix detectors

The Timepix detector used for this project is the member of the 23 Medipix/Timepix detector family^[32-35]. These are active pixels detectors 24 developed by the Medipix collaboration at the European Organization for 25 Nuclear Research (CERN, Geneva, Switzerland). They are based on 26 Complementary Metal-Oxide Semiconductor (CMOS) technology^[36]. The 27 Timepix application-specific integrated circuit (ASIC) is an improved 28 version of the Medipix ASIC with two new measurement modes in addition 29 to the simple particle counting mode. One is the Time-over-Threshold 30 mode (TOT mode) in which every pixel registers how long it is above a 31 certain charge threshold level. The other operating mode, and the one 32 used for mass spectrometry, is the Time-of-Flight mode (TOF mode). One 33 Timepix chip consists of 256×256 pixels, each capable of individually 34 measuring the time-of-flight of ions hitting the detector with respect to an 35 externally applied trigger signal. Each clock cycle is 10 ns wide and the 36 maximum measurement interval of the Timepix detector is 118 µs. The 37 size of an individual pixel is 55 µm×55 µm. It is also possible to build 2×2n 38 39 arrays of Timepix chips if bigger detectors are necessary for the experiments. 40

41 Microscope mode mass spectrometer

A TRIple Focusing Time-of-Flight (TRIFT II) mass spectrometer 1 (Physical Electronics, Inc., Chanhassen, MN, USA) is used. The mass 2 spectrometer is equipped with a 20 keV C₆₀ primary ion source (lonoptika, 3 Chandlers Ford, Hampshire, United Kingdom) and a high voltage Timepix 4 detector setup^[29]. The primary ion source is operated at 20 keV primary 5 ion energies with the C_{60}^{2+} ion selected for high primary ion beam current. 6 The primary ion beam uses a pulse length of 60 ns and is then bunched 7 for better spectral resolution. All the measurements are done in static 8 SIMS mode, where the primary ion dose is well below the static limit 9 $(1x10^{13} \text{ ions/cm}^2).$ 10

The high voltage Timepix setup installed on the instrument has 11 been described in details elsewhere^[29]. Briefly, it consists of a 2×2 array of 12 Timepix chips behind a chevron MCP stack. The MCPs are operated at a 13 bias of 1.375 kV for the positive mode and 1.4 kV for the negative mode 14 experiments, unless it is stated otherwise. The data readout system uses 15 the ReLAXD (high Resolution Large Area X-ray Detector) readout board 16 which has a readout speed of 1 Gbit/s^[37, 38] and is operated at a frame rate 17 of 10 frames/s. The chips are cooled with a Peltier element based active 18 cooling system. The entire detector setup can be floated at +12 kV (for 19 negative ion mode) or -8 kV (for positive ion mode) with the use of the 20 TRIFT II mass spectrometer's high voltage power supply. This offers 21 additional post-acceleration of the secondary ions before they reach the 22 MCPs. Both the Timepix/RelaXD system and the C₆₀ source are triggered 23 from the TRIFT II mass spectrometer master trigger, which is down 24 sampled 100 times to 10 Hz. For the data acquisition the Pixelman data 25 acquisition software is used^[39]. 26

27 The Timepix detector is operated in TOF mode. In all modes, a typical, sparse data frame contains the x- and y-coordinate of every 28 triggered pixel. Additionally, in TOF mode, the data file lists the time-of-29 flight. In the TOF mode, every measurement frame contains the TOF 30 information obtained from a single primary ion pulse. The mass spectrum 31 is reconstructed by making a histogram of the TOF values from the 32 separate frames. Standards are used to calculate mass calibration 33 parameters, which are then applied to the TOF spectra. Total ion images 34 are constructed by summing all of the individual frames. Selected ion 35 images are plotted by extracting the pixel positions and intensities for a 36 selected mass spectral peak. 37

38


39 Samples

40 Various benchmark samples are used to test the imaging and mass 41 spectral performance of the system. These include brilliant green dye

(green Staedtler Lumocolor 318-5 permanent marker, Staedtler Mars 1 GmbH & Co. KG, Nuernberg, Germany) as well as a 1 mg/mL solution of a 2 mixture of different chain length polyethyleneglycols (PEG 200-3500) 3 mixed with 7 mg/mL α -cyano-4-hydroxycinnamic acid (CHCA) in 1:1 ratio 4 (1 µL deposited). The benchmark samples are placed on an Indium-Tin 5 6 oxide (ITO) coated glass slide (4-8 Ω resistance, Delta Technologies, 7 Stillwater, MN, USA). A hexagonal transmission electron microscopy (TEM) grid (700 mesh, G2760N, 3.05 mm diameter, 37 µm pitch, 8 µm bar 8 width; Agar Scientific Limited, Stansted, United Kingdom) is placed on the 9 top of the samples. For biological tissue imaging, 12 um thick coronal 10 mouse brain (male balb/c mouse; Harlan Laboratories, Boxmeer, The 11 Netherlands) sections are used. The brain is sectioned in a Microm 12 HM525 cryomicrotome (Thermo Fisher Scientific, Walldorf, Germany) and 13 the sections are placed on an ITO coated glass slide. The sections are 14 15 kept at -20 °C until further use. Before measurement the samples are dried in a vacuum desiccator and are measured without any further sample 16 preparations steps. 17

18

19 **Results and discussion**

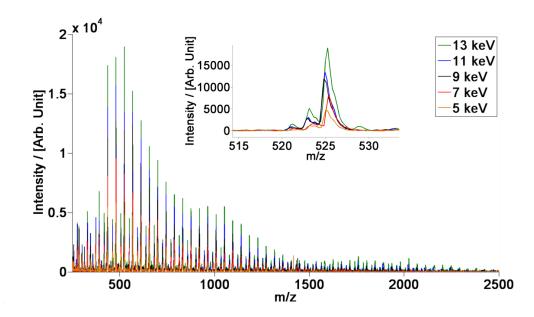

20

Figure 1. Mass spectra of a mixture of polyethylene glycol measured with C_{60} SIMS on the Timepix detector at different ion acceleration energies. The inset shows a selected PEG peak at m/z 525. Higher ion acceleration energies improve the signal-to-noise ratio and extend the accessible upper mass range.

25

As a first experiment the system's spectral quality is assessed. Also, the effect of the increased post-acceleration of the ions on the mass spectra is

systematically studied. A mixture of PEGs is used and 10,000 frames at 1 different total ion acceleration energy values between 5 keV and 13 keV. 2 The MCP is operated at an MCP bias of 1.5 kV. Figure 1. shows an 3 overlay of the PEG mass spectra obtained at the different total ion 4 acceleration energies. The sodiated molecular ions of the PEGs are 5 6 detected in this measurement. Increasing the acceleration energy of the secondary ions has a two-fold effect on the mass spectra. First, the 7 intensity of the detected ions increases. Thus, the higher post acceleration 8 yields a better signal-to-noise ratio (S/N of 184 for the ion at m/z 525 at 13 9 keV total acceleration energy and S/N 55 at 5 keV total acceleration 10 11 energy). Also, the accessible mass range is significantly increased at the highest energy values. This is due to the higher ion energies, where more 12 ions have a chance to start an electron cascade in the MCP. Also, higher 13 mass ions that would not have the necessary energy to start an electron 14 15 cascade can be detected due to the additional post-acceleration capabilities of the system. It is important to note that the same sample 16 measured for 10,000 frames with an Au primary ion source with a higher 17 MCP gain results in a spectrum with a S/N of 37, which is five times lower 18 19 than the S/N achieved with the combination of the C_{60} source and the high voltage Timepix setup, and mass range between m/z 0 and 1200 (see 20 supplementary figure S1). Thus, the use of the C_{60} primary ion source 21 results in a significant increase in spectral guality, even at low post-22 23 acceleration, due to the higher secondary ion yields associated with such 24 polyatomic primary ion sources.

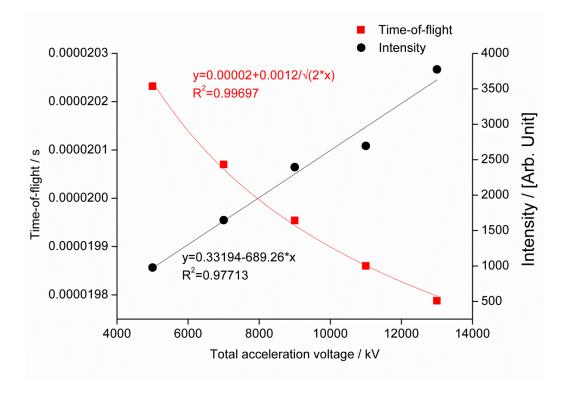


Figure 2. Intensity (black) and the time-of-flight as the number of clock cycles (red) of the PEG ion at *m/z* 526 as a function of the ion acceleration energy.

The high voltage setup has two main effects on the spectral quality. 4 The first one is the increase of the secondary ion intensities discussed in 5 6 the previous paragraph. The second effect is a decreasing time of flight of the same ion as a function of the acceleration voltage. To demonstrate 7 these effects the measurements shown in Fig. 1 are divided into 5 x 2,000 8 9 frame segments. Figure 2. shows the average intensities and average 10 time-of-flight values of these five smaller datasets as a function of the ion acceleration energies. The time-of-flight of the selected ion changes as a 11 quadratic function of the acceleration voltage. This is in agreement with 12 the TOF analyzers calibration equation (1) 13

$$t = \frac{d}{\sqrt{2U}} \sqrt{\frac{m}{q}}$$

1

14

(1)

where t is the time-of-flight of the ion, d is the length of the ion's flight path, m is the mass of the ion, q is the charge of the ion and U is the acceleration voltage. Analysis of the dependence of ion intensity demonstrates a linear increase with ion kinetic energy, as expected.

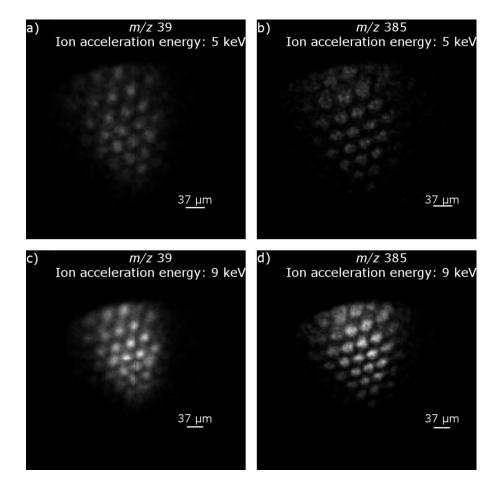




Figure 3. Comparison of the image quality at different ion acceleration energies. Selected-ion images of potassium at *m/z* 39 (a, c) and brilliant green dye at *m/z* 385 (b, d) are plotted on the same intensity scale. The detector is operated at total ion acceleration energies of 5 keV (a, b) or 9 keV (c, d).

6 Figure 3 demonstrates the imaging capabilities of the system on 7 brilliant green dye underneath a TEM grid, in the positive ion mode. The same area of the sample is measured with the detector at ground (Figure 8 3a and 3b) and at 4 kV (Figure 3c and 3d). The selected-ion images for 9 the comparison are the cation of the brilliant green dye (m/z 385) and 10 potassium ion (m/z 39) which both localize in the holes of the grid. The 11 size of the primary ion beam is \sim 224 µm × 272 µm and the images are 12 reconstructed from 20,000 frames. Fig. 3 shows that the high voltage 13 setup is capable to provide the same image quality as the previous 14 detector setup (where the detector is held at ground potential). The effects 15 of the higher ion acceleration on the image guality are a better image 16 contrast, due to the higher S/N ratio, and a smaller magnification factor, 17 due to the shorter time that the ions spend in the magnification region of 18 the instrument. The observed spatial resolving power is ~7 µm (see 19 supplementary figure S2.) with a pixel size of 900 nm. This value is close 20 to the previously reported^[31] spatial resolutions for the Timepix setup 21 without high voltage capabilities and for a microprobe mode reference 22

- 1 measurement system that uses the combination of an MCP with a time-to-
- 2 digital converter (TDC), both on a TRIFT mass spectrometer.

3

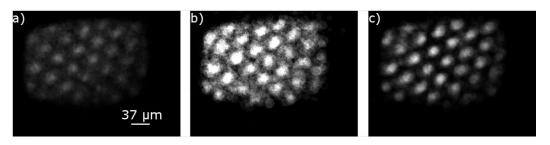


Figure 4. Negative mode imaging of brilliant green dye under a TEM grid. Selected-ion images show the distribution of C₂H at m/z 25 (a), C₄H at m/z 49 (b) and HSO₄ at m/z 97 (c).

The previous Timepix setup was only operable in positive ion mode. 7 One of the advantages of the new high voltage setup is the capability to 8 9 detect negative ions as well. Thus, it is capable of providing complementary information to the positive mode measurements. This new 10 functionality of the system was previously demonstrated on peptide and 11 protein standards with MALDI ionization^[29]. However, no negative mode 12 images were recorded thus far with a pixelated detector. Figure 4. shows 13 the first example of negative mode microscope mode SIMS imaging with a 14 Timepix detector. The sample is the same grid standard used for the 15 positive mode measurements shown in Fig. 3. The detector is operated at 16 an ion acceleration energy of 12.5 keV. As expected, the ions detected are 17 mostly the standard small negative ions common to SIMS, such as $[C_2H]^{-1}$ 18 at m/z 25, $(C_4H)^-$ at m/z 49; or $(PO_3)^-$ at m/z 79. These ions are abundant 19 in the holes of the sample which means they are present in the green dye 20 underneath the grid. The image quality is similar to the positive mode 21 22 image quality as can be seen in supplementary Figure S2. This evaluation reveals a resolving power of 7.5 µm in negative mode which is in 23 accordance with the images acquired in positive ion mode. 24

The fivefold increase in sensitivity achieved with the new setup 25 makes the imaging of biological samples possible. Figure 5 exemplifies 26 27 the system's improvement for the imaging of biological samples. This first 28 example of a mosaic microscope mode SIMS imaging with a pixelated detector shows complementary positive and negative ion mode images of 29 a biological tissue sample. In particular, the sample is half of the anterior 30 commissure area of a mouse brain section as can be seen on the 31 32 hematoxylin and eosin (H&E) stained image (Figure 5a). The images are reconstructed from 4 tiles. The step size between each tile is 100 µm, 33 such that there is sufficient overlap between the tiles to correct for the 34 35 elliptic shape of the primary ion beam spot during the reconstruction of the 36 image. In both positive and negative mode, 20,000 frames per tile are 37 collected. One measurement takes 80,000 frames compared to the single 10

tile imaging with the gold primary ion source that took more than 300,000 1 2 frames to achieve similar image quality. This means that as a result of the 3 improvements, ~80% fewer frames are needed to measure an area that is roughly four times the size as the previously reported microscope mode 4 SIMS image with the Timepix detector. Further, the increased secondary 5 6 ion yield of the C_{60} source reduces the measurement time to only half of 7 the time reported earlier for the single tile microscope mode SIMS measurement. Also, intact phospholipids are detected from the tissue (see 8 supplementary figure S3). 9

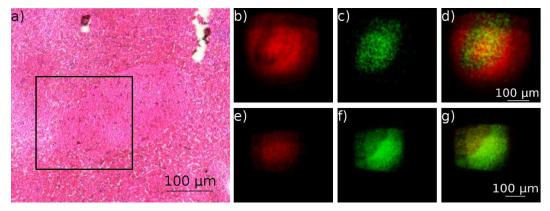


Figure 5. Microscope image of a H&E stained mouse brain (a), positive (b, c, d) and negative (e, f, g) microscope mode SIMS image of the anterior commissure region of the mouse brain. The selected ion images show the distribution of sodium at m/z 23 (b), the cholesterol fragment $[M+H-H_2O]^+$ at m/z 369 (c), CH⁻ at m/z 13 (e) and PO₃⁻ at m/z 79 (f) and the overlay of the two positive (d) and negative mode (g) images

It is possible to distinguish the anterior commissure in both positive 17 and negative mode. The positive mode selected ion images show the 18 distribution of the sodium ion at m/z 23, which has higher abundance in 19 20 the tissue surrounding the anterior commissure and a cholesterol fragment 21 $([M+H-H_2O]^+)$ at m/z 369, which localizes in the anterior commissure. The ions selected for the negative mode images are $(CH)^{-}$ at m/z 13 with a 22 23 homogenous distribution in the imaged area of the tissue and (PO₃) at m/z79, which is localized in the tissue around the anterior commissure. 24

25

10

26 Conclusions

For the first time, the combination of a C_{60} primary ion source and a pixelated detector system with ion post acceleration capabilities was used for high spatial resolution microscope mode SIMS imaging in both positive and negative ion mode. The combination of the higher secondary ion yield of the C_{60} primary ion source and the higher ion post acceleration possibilities of the new Timepix setup resulted in an increased signal-tonoise ratio and wider mass range compared to the earlier microscope mode SIMS experiments. The acquisition time necessary to achieve high quality data was also significantly reduced. The image quality, namely the spatial resolving power, is comparable to the earlier Timepix based microscope mode SIMS studies. It is also possible, for the first time, to measure microscope mode SIMS images in negative mode with a pixelated detector.

Timepix based microscope mode SIMS imaging opens up new 8 9 possibilities, such as the fast 3D imaging of complex samples with subcellular spatial resolution. However, several challenges such as the speed 10 of the readout system and the single stop nature of the pixels remain that 11 12 require further improvements on the detector level. These include a faster, 1 kHz readout system, multi-hit capabilities and a new, compact data 13 format. Also, online, on-the-fly data analysis capability integrated on the 14 chips or on-board the readout electronics would be advantageous for 15 16 future applications and are subject to present and future studies. Microscope mode SIMS imaging with a pixelated detector shows 17 promising results for high resolution mass spectrometry imaging of 18 biological systems. 19

20

21 ACKNOWLEDGEMENTS

This work is part of the research program of the Foundation for 22 Fundamental Research on Matter (FOM), which is part of the Netherlands 23 Organisation for Scientific Research (NWO). Part of this research is 24 25 supported by the Dutch Technology Foundation STW, which is the applied science division of NWO, and the Technology Program of the Ministry of 26 Economic Affairs, Project OTP 11956. This publication was supported by 27 the Dutch national program COMMIT and the Netherlands Proteomics 28 29 Center. The authors acknowledge Sheng Chen (Physical Electronics) for excellent support during this project by sharing his knowledge on the 30 TRIFT II system and Paul Blenkinsopp and Rowland Hill (Ionoptika) for 31 their help with C₆₀ primary ion source. The authors are thankful to Berta 32 33 Cillero Pastor who performed the pathological staining of the tissue samples, Frans Giskes and Ronald Buijs for their technical support and 34 Ivo Klinkert for his support with the data processing. 35

1 References

- 2 [1] K. Chughtai, R. M. A. Heeren, *Chemical Reviews* **2010**, *110*, 3237.
- 3 [2] L. A. McDonnell, R. M. A. Heeren, *Mass Spectrometry Reviews* 2007, 26, 606.
- 4 [3] T. L. Colliver, C. L. Brummel, M. L. Pacholski, F. D. Swanek, A. G. Ewing, N. 5 Winograd, *Anal Chem* **1997**, *69*, 2225.
- 6 [4] A. F. M. Altelaar, I. Klinkert, K. Jalink, R. P. J. de Lange, R. A. H. Adan, R. M. A.
 7 Heeren, S. R. Piersma, *Anal Chem* 2006, *78*, 734.
- 8 [5] B. Cillero-Pastor, G. Eijkel, A. Kiss, F. J. Blanco, R. M. A. Heeren, *Anal Chem* 2012,
 9 84, 8909.
- [6] D. Debois, K. Hamze, V. Guerineau, J. P. Le Caer, I. B. Holland, P. Lopes, J.
 Ouazzani, S. J. Seror, A. Brunelle, O. Laprevote, *Proteomics* 2008, *8*, 3682.
- 12 [7] N. Desbenoit, I. Schmitz-Afonso, C. Baudouin, O. Laprevote, D. Touboul, F.
 13 Brignole-Baudouin, A. Brunelle, *Analytical and bioanalytical chemistry* 2013, 405, 4039.
- [8] L. Fornai, A. Angelini, I. Klinkert, F. Giskes, A. Kiss, G. Eijkel, E. A. Amstalden-van
 Hove, L. A. Klerk, M. Fedrigo, G. Pieraccini, G. Moneti, M. Valente, G. Thiene, R.
 M. A. Heeren, *Analytical and bioanalytical chemistry* 2012, 404, 2927.
- 18 [9] L. A. Klerk, P. Y. W. Dankers, E. R. Popa, A. W. Bosman, M. E. Sanders, K. A.
 19 Reedquist, R. M. A. Heeren, *Anal Chem* **2010**, *82*, 4337.
- 20 [10] H. Liebl, Journal of Applied Physics **1967**, *38*, 5277.
- 21 [11] S. Chandra, D. R. Smith, G. H. Morrison, *Anal Chem* **2000**, *72*, 104 A.
- I. V. Veryovkin, W. F. Calaway, C. E. Tripa, M. J. Pellin, Nuclear Instruments &
 Methods in Physics Research Section B-Beam Interactions with Materials and
 Atoms 2007, 261, 508.
- 25 [13] N. Winograd, Anal Chem **2005**, 77, 142a.
- [14] D. Weibel, S. Wong, N. Lockyer, P. Blenkinsopp, R. Hill, J. C. Vickerman, *Anal Chem* 2003, *75*, 1754.
- [15] J. Xu, S. Ostrowski, C. Szakal, A. G. Ewing, N. Winograd, *Applied Surface Science* **2004**, 231, 159.
- 30 [16] R. Kersting, B. Hagenhoff, F. Kollmer, R. Mollers, E. Niehuis, *Applied Surface*31 Science 2004, 231, 261.
- 32 [17] A. M. Piwowar, J. C. Vickerman, *Surface and Interface Analysis* **2010**, *42*, 1387.
- 33 [18] D. E. Weibel, N. Lockyer, J. C. Vickerman, *Applied Surface Science* **2004**, *231*, 146.
- S. C. C. Wong, R. Hill, P. Blenkinsopp, N. P. Lockyer, D. E. Weibel, J. C. Vickerman,
 Applied Surface Science 2003, 203, 219.
- 36 [20] A. Wucher, *Applied Surface Science* **2006**, *252*, 6482.
- S. L. Luxembourg, L. A. McDonnell, T. H. Mize, R. M. A. Heeren, *Journal of Proteome Research* 2005, *4*, 671.
- 39 [22] S. L. Luxembourg, T. H. Mize, L. A. McDonnell, R. M. A. Heeren, *Anal Chem* 2004,
 40 *76*, 5339.
- 41 [23] L. A. Klerk, N. P. Lockyer, A. Kharchenko, L. MacAleese, P. Y. W. Dankers, J. C.
 42 Vickerman, R. M. A. Heeren, *Anal Chem* 2010, *82*, 801.
- 43 [24] M. Froesch, S. L. Luxembourg, D. Verheijde, R. M. A. Heeren, *Eur J Mass Spectrom* 2010, *16*, 35.
- 45 [25] C. Bamberger, U. Renz, A. Bamberger, J Am Soc Mass Spectr 2011, 22, 1079.
- 46 [26] J. H. Jungmann, A. Gijsbertsen, J. Visser, J. Visschers, R. M. A. Heeren, M. J. J.
 47 Vrakking, *Rev Sci Instrum* **2010**, *81*.
- 48 [27] J. H. Jungmann, L. MacAleese, R. Buijs, F. Giskes, A. de Snaijer, J. Visser, J.
 49 Visschers, M. J. J. Vrakking, R. M. A. Heeren, J Am Soc Mass Spectr 2010, 21,
 50 2023.

- [28] J. H. Jungmann, L. MacAleese, J. Visser, M. J. J. Vrakking, R. M. A. Heeren, *Anal Chem* 2011, *83*, 7888.
- J. H. Jungmann, D. F. Smith, A. Kiss, L. MacAleese, R. Buijs, R. M. A. Heeren,
 International Journal of Mass Spectrometry 2013, 341–342, 34.
- [30] J. H. Jungmann, D. F. Smith, L. MacAleese, I. Klinkert, J. Visser, R. M. A. Heeren, J
 Am Soc Mass Spectr 2012, 23, 1679.
- 7 [31] A. Kiss, J. H. Jungmann, D. F. Smith, R. M. A. Heeren, *Rev Sci Instrum* **2013**, *84*.
- 8 [32] X. Llopart, R. Ballabriga, M. Campbell, L. Tlustos, W. Wong, Nuclear Instruments
 9 & Methods in Physics Research Section a-Accelerators Spectrometers Detectors
- 10 *and Associated Equipment* **2007**, *581*, 485.
- [33] X. Llopart, M. Campbell, Nuclear Instruments & Methods in Physics Research
 Section a-Accelerators Spectrometers Detectors and Associated Equipment 2003,
 509, 157.
- 14 [34] X. Llopart, M. Campbell, R. Dinapoli, D. S. Segundo, E. Pemigotti, *leee Transactions on Nuclear Science* 2002, *49*, 2279.
- [35] M. Campbell, M. Collaboration, Nuclear Instruments & Methods in Physics
 Research Section a-Accelerators Spectrometers Detectors and Associated
 Equipment 2011, 633, S1.
- 19 [36] F. M. Wanlass, US patent # 3,356,858, **1967**.
- [37] J. Visser, B. van der Heijden, S. J. A. Weijers, R. de Vries, J. L. Visschers, Nuclear
 Instruments & Methods in Physics Research Section a-Accelerators
 Spectrometers Detectors and Associated Equipment 2011, 633, S22.
- [38] Z. Vykydal, J. Visschers, D. S. Tezcan, K. De Munck, T. Borgers, W. Ruythooren, P.
 De Moor, Nuclear Instruments & Methods in Physics Research Section a Accelerators Spectrometers Detectors and Associated Equipment 2008, 591, 241.
- [39] T. Holy, J. Jakubek, S. Pospisil, J. Uher, D. Vavrik, Z. Vykydal, Nuclear Instruments
 & Methods in Physics Research Section a-Accelerators Spectrometers Detectors
 and Associated Equipment 2006, 563, 254.
- 29
- 30
- 31
- 32
- 33

1 Cluster SIMS Microscope Mode Mass Spectrometry 2 Imaging

- 3 András Kiss¹, Donald F. Smith¹, Julia H. Jungmann¹, Ron M.A. Heeren^{1*}
- ¹ FOM Institute AMOLF, Science Park 104, 1098 XG Amsterdam, The
 5 Netherlands


6 * Author to whom correspondence should be addressed. Email: heeren@amolf.nl

7


8

9

SUPPLEMENTARY FIGURES

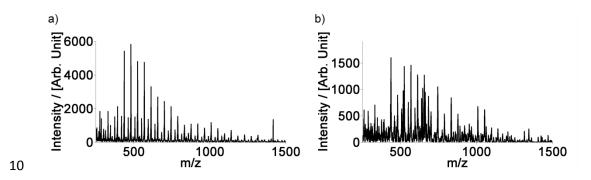


Supplementary figure S1. Comparison of a PEG spectra recorded with the C_{60}^{+} Timepix setup (a) and with the Au⁺ LMIG Timepix setup (b)

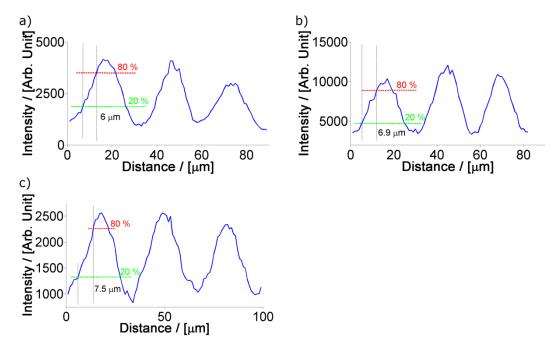

¹³

Supplementary figure S2. Line scans from the total ion image of the TEM grid in positive mode at a total acceleration energy of 5 keV (a) and 9 keV (b) and in negative mode at the total acceleration energy of 12.5 keV (c)

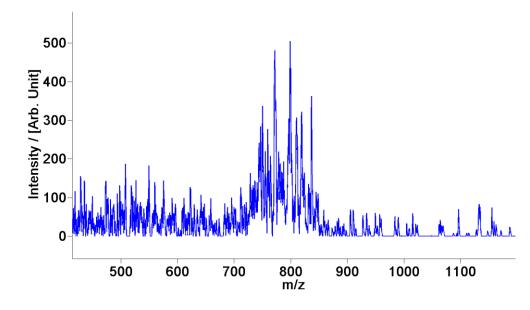

Supplementary figure S3. Mass spectrum from the positive mode tissue imaging experiment. The
 spectrum shows the detected intact phospholipids from the tissue