

# Fast and accurate alignment of long bisulfite-seq reads

Brent S. Pedersen<sup>1,\*</sup>, Kenneth Eyring<sup>1</sup>, Subhajyoti De<sup>1,2</sup>, Ivana V. Yang<sup>1</sup> and David A. Schwartz<sup>1</sup>

<sup>1</sup>Department of Medicine, University of Colorado Denver, School of Medicine, Denver, Colorado, USA. 80045

<sup>2</sup>University of Colorado Cancer Center, Molecular Oncology Program, Aurora, Colorado, USA

Received on XXXXX; revised on XXXXX; accepted on XXXXX

Associate Editor: XXXXXXXX

## ABSTRACT

**Summary:** We introduce a new tool, *bwa-meth*, to align bisulfite-treated sequences and compare it to existing aligners. We show that it is fast and accurate even without quality-trimming and the output is immediately usable by downstream tools. We gauge accuracy by comparing the number of on and off-target reads from a targeted sequencing project and by simulations.

**Availability and Implementation:** The benchmarking scripts and the *bwa-meth* software are available at <https://github.com/brentp/bwa-meth/> under the MIT License.

**Contact:** bpederse@gmail.com

**Supplementary information:** Supplemental Information 1

## 1 INTRODUCTION

Bisulfite sequencing (BS-Seq) is a common way to explore DNA methylation. As a result, software has been developed to map sequence reads treated with bisulfite to a reference genome (Frith *et al.*, 2012; Pedersen *et al.*, 2011; Wu and Nacu, 2010; Krueger and Andrews, 2011; Xi and Li, 2009; Hansen *et al.*, 2012). Previous studies have compared alignment statistics on real (Pedersen *et al.*, 2011; Xi and Li, 2009; Shrestha and Frith, 2013) and simulated (Frith *et al.*, 2012) reads, however most of these are limited by knowledge of the ground-truth and assumptions of the simulation, respectively.

Here, we present an analysis of current BS-Seq mappers including the "four-base" aligners Last (Frith *et al.*, 2012), GSNAP (Wu and Nacu, 2010) and BSMAP (Xi and Li, 2009) as well as the "three-base" aligners BSmooth (Hansen *et al.*, 2012), Bison, and Bismark (Krueger and Andrews, 2011) which perform *in silico* conversion of cytosines to thymines. In addition, we introduce our own, simple three-base aligner that wraps *BWA mem* (Li, 2013).

The comparison is performed on 100-base paired-end reads which are of modest length by current standards, but, to our knowledge, longer than utilized in any comparison. We hypothesized that long, paired reads, with up to 200 bases from the same genomic region, could change the decision on which alignment method performed the best and that aligners which relied on global alignment might have reduced performance on longer reads. We found limitations

to existing aligners including the writing of large temporary files, required quality-trimming, high memory-use, long run-time, output that was not suitable for consumption by traditional tools, or some combination of these inconveniences. We wrote a BS-Seq aligner based on BWA mem (Li, 2013) to address these limitations. This new aligner, *bwa-meth*, allows indels, local alignments, and it never writes a temporary-file of the reads to disk, instead streaming the *in silico* converted reads directly to the aligner and streaming the alignments directly to a well-formatted alignment file suitable for use in downstream tools. In addition, it works well without quality-trimming, thereby reducing storage requirements 3-fold by not creating quality-trimmed or converted reads.

## 2 APPROACH


In order to determine the accuracy of an aligner, we utilize a dataset from Agilent's SureSelect Mouse Methyl-Seq kit which captures about 99 million bases from CpG dense regions in the mouse genome (a similar approach is available for human regions). We evaluate an aligner by the number of reads in the capture area as compared to outside of the capture area. While there will be off-target capture, all aligners are subject to the same assumptions. With those constraints, we can plot a receiver operating curve (ROC) with true positives as reads within and false positives as reads outside of the target regions. Shrestha and Frith, 2013 performed an unbiased analysis on 35 base reads sequenced from the X chromosome and called a read as accurately mapped if it aligned to the X chromosome; here, we map to the entire genome with 100-base paired-end reads.

In addition, we map 100-base paired-end indel-free reads simulated using Sherman (v0.1.6), the software from the authors of Bismark. All data were aligned to mouse genome version *mm10*.

## 3 METHODS

We aligned real and simulated data, both as-is and trimmed by quality (using Trim Galore [[http://www.bioinformatics.babraham.ac.uk/projects/trim\\_galore/](http://www.bioinformatics.babraham.ac.uk/projects/trim_galore/)] default parameters), using the software and versions in Supplemental Table 1. We evaluated a few parameters for each method and report only the best-performing here. Likely, every aligner could show improved results with exhaustive search of the parameters, but this is a representation of reasonable

\*To whom correspondence should be addressed. Email: bpederse@gmail.com

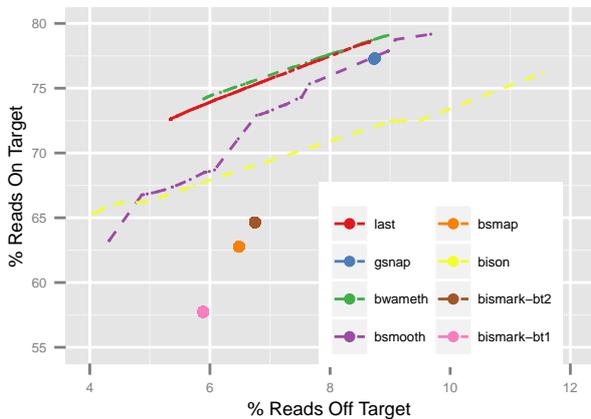

**Fig. 1.** Percent of un-trimmed paired-end, 100-base reads on (y) and off (x) target for the tested aligners. Aligners that report mapping quality are shown as connected dots for each quality cut-off. Reads are limited to those considered as primary, mapped alignments by the aligner. Note that, as shown in the supplement, many aligners have better performance on trimmed reads.

selections. We considered a real read to be on-target if it was within 1001 bases of a target region.

We designed *bwa-meth* for paired-end reads from the directional protocol but it can align single-end reads. *bwa-meth* outputs alignments directly to a BAM file that is usable even by those tools that require coordinate-sorted alignments and read-groups. Since it consists of fewer than 600 lines of code and runs quickly, it can be used as a platform to test other optimizations. For example in Supplemental Figure 7, we show a feature that further reduces the number of off-target reads by considering only the single strand targeted by the SureSelect protocol.

Although our comparisons are on an inbred mouse strain with few polymorphisms we expect the results will hold even with insertions and deletions due to our use of BWA mem (Li, 2013).

## 4 DISCUSSION

### 4.1 Accuracy

For the aligners that report a range of mapping quality scores—an indicator of the aligner’s confidence in the alignment—we vary the score from 1 to the maximum, 255, to draw an ROC-like curve showing the trade-off between sensitivity and specificity. For the other aligners, we plot their single location. Figure 1 shows the on and off-target reads for our real paired-end data. Last and *bwa-meth* align the most reads on target with a low percent of off-target reads, but Last can provide better control over the number of off-target reads. On trimmed data, *bwa-meth*, Last, Bison and Bsmooth all perform similarly with Bison allowing very good control over the number off-target reads (Supplement Figure 2). Only Last, GSNAP, and *bwa-meth* are relatively unaffected by quality-trimming.

Bismark performance varies depending on whether it uses bowtie version 1 or 2 for the backend but in either case, (as with all of the bowtie-backed aligners) it performs better with trimmed reads.

A similar comparison is shown in Supplemental Figure 2 for trimmed data. We also show the ROC figures for simulated data with errors in Supplemental Figures 3 (original) and 4 (trimmed)

and without errors in Supplemental Figures 5 and 6. *Bwa-meth* does out-perform all aligners for simulated data.

### 4.2 Computational Resources

We are more interested in the accuracy of a method than the speed. In general the aligners are comparable in terms of speed. While Bismark with bowtie1 is quite fast on the simulated data it can not use multiple processes and so takes the longest user-time. Bsmmap takes the longest CPU time on the real data. We report the exact timings and maximum memory use in Supplemental Information. Bsmmap uses the least disk, never writing an index of the reference genome and only writing the alignment files. All other aligners write an index of the reference genome. Last, bsmooth, and bismark write additional copies of the reads to disk. *Bwa-meth* avoids writing the *in silico* converted reads to disk by streaming them directly to the aligner and shows nearly identical accuracy without read-trimming. The trimming and conversion steps can increase data storage needs enough to be a burden in our experience.

None of the programs used an inordinate amount of memory, however, due to the parallelization strategy, Last did require about 10GB of shared memory per process.

## 5 CONCLUSION

We have used simulated reads along with reads from a capture method that targets CpG-rich regions to compare aligners. We show that *bwa-meth* is very accurate, even without quality trimming. Aligners that use global alignment benefit more from quality trimming and may also require an additional copy of the *in silico* converted reads. *bwa-meth* is also fast and it outputs alignments to a sorted bam that is immediately usable by downstream tools.

## ACKNOWLEDGEMENT

We thank Devon Ryan for several helpful conversations.

**Funding:** This was funded by R01 HL097163, R01 HL101251, 1I01BX001534, RC2 HL101715, N01 AI90052, and S10 RR031832.

## REFERENCES

- Frith, M. C., Mori, R., and Asai, K. (2012). A mostly traditional approach improves alignment of bisulfite-converted DNA. *Nucleic acids research*, **40**(13), e100–e100.
- Hansen, K. D., Langmead, B., Irizarry, R. A., et al. (2012). Bsmooth: from whole genome bisulfite sequencing reads to differentially methylated regions. *Genome Biol*, **13**(10), R83.
- Krueger, F. and Andrews, S. R. (2011). Bismark: a flexible aligner and methylation caller for bisulfite-seq applications. *Bioinformatics*, **27**(11), 1571–1572.
- Li, H. (2013). Aligning sequence reads, clone sequences and assembly contigs with bwa-mem. *arXiv preprint arXiv:1303.3997*.
- Pedersen, B. S., Hsieh, T.-F., Ibarra, C., and Fischer, R. L. (2011). MethylCoder: software pipeline for bisulfite-treated sequences. *Bioinformatics*, **27**(17), 2435–2436.
- Shrestha, A. M. S. and Frith, M. C. (2013). An approximate bayesian approach for mapping paired-end dna reads to a reference genome. *Bioinformatics*, **29**(8), 965–972.
- Wu, T. D. and Nacu, S. (2010). Fast and SNP-tolerant detection of complex variants and splicing in short reads. *Bioinformatics*, **26**(7), 873–881.
- Xi, Y. and Li, W. (2009). BSMAP: whole genome bisulfite sequence MAPPING program. *BMC bioinformatics*, **10**(1), 232.

# Supplemental Information

Pedersen *et al.*

## 1 Advantages

*bwa-meth* provides a number of advantages:

1. ease of use – a single script indexes the fasta, aligns the reads, and tabulates methylation. Default parameters work well.
2. speed – due to the efficient parallelization by *bwa*, *bwa-meth* can run on as many CPU's are available on a single machine.
3. disk-usage – *bwa-meth* can take either compressed or uncompressed files and the reads are streamed directly to the aligner, never written to disk. An aligner that requires trimmed reads and that writes the converted reads to disk will require 3X as much storage for just the raw sequence data.
4. useful output – the BAM that results from *bwa-meth* is sorted, indexed and contains the proper mapping quality scores, alignment flags, and read-groups. In addition, the header @SQ lines are sorted so that tools like Picard and GATK will accept them as is.
5. simplicity – we have carefully designed *bwa-meth* to be quite simple and the code readable so that enhancements are simple. For example, see #6.
6. strand-specificity – Devon Ryan, the author of Bison, noted that many capture methods target only the reverse strand – commonly called the original bottom. Though this would be simple to add to any aligner, ours is the only one that currently supports this. See from the supplemental figure 5 for the reduction in off-target reads that results from only considering reads that align to the original bottom strand.
7. specificity and sensitivity – while the difference between aligners in this regard is small, we have shown that *bwa-meth* consistently performs as the best or among the best. We believe this to be the case in general and not limited to the specific cases that we have shown.
8. accuracy without trimming – the accuracy is due in part to *bwa*'s local alignment. This is something that is not supported in *bismark* and means that we can achieve very good accuracy without performing additional trimming which would result in additional disk-usage and processing time.

## 2 Aligners Compared

We used a modified version of the calling script for last available in the github repository.

Table 1: Methods Compared

software	version	arguments
bismark (bowtie1)	0.11.1	-gzip -bam -n 2 -l 24
bismark (bowtie2)	0.11.1	-gzip -bam -N 1 -score_min L,-0.4,-0.5
bsmap	2.74	bsmap -v3 -m0 -x1000 -S42 -n0 -s12 -I1
bsmooth	fb3f7ef	-very-sensitive
bison	0.3.0	-directional -very-sensitive-local -N 1
bwa-meth	0.10	bwa-meth
gsnap	2014-02-28	-npaths 1 -quiet-if-excessive -k 15
last	392	last-bisulfite-paired.sh

### 3 Comparison On Real Data

We compared 7 aligners on real data as described in the text. Supplemental Figure 1 below shows the same data as Figure 1 from the main text.

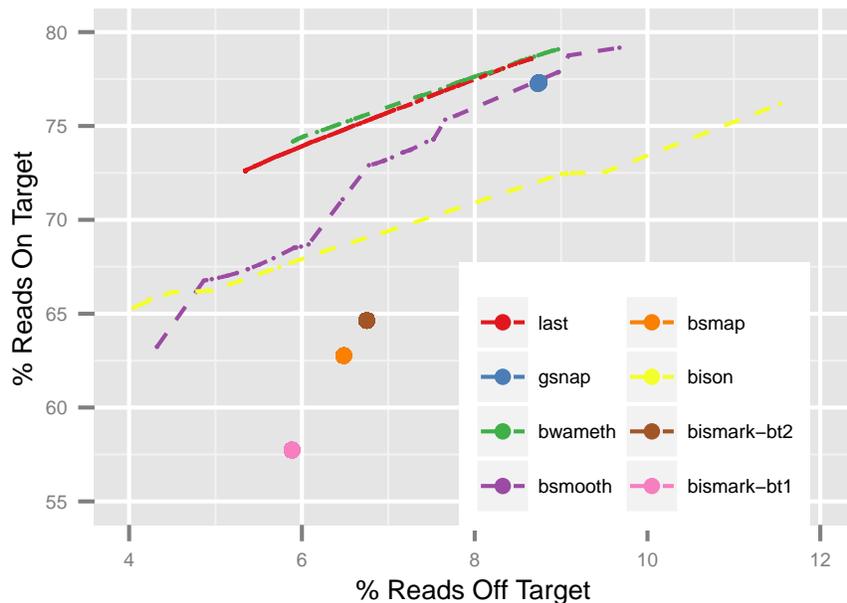


Figure 1: Percent of paired-end, 100-base reads on (y) and off (x) target for the tested aligners. Aligners that report mapping quality are shown as lines that span each quality cut-off. Reads are limited to those considered as primary, mapped alignments by the aligner. This is a color version of figure 1 from the paper

When we trim the reads with *trim\_galore* before aligning, the result is shown below in Supplemental Figure 2. We report the percent of reads aligned relative to the original count in the untrimmed data because we are interested in the overall mapping rate.

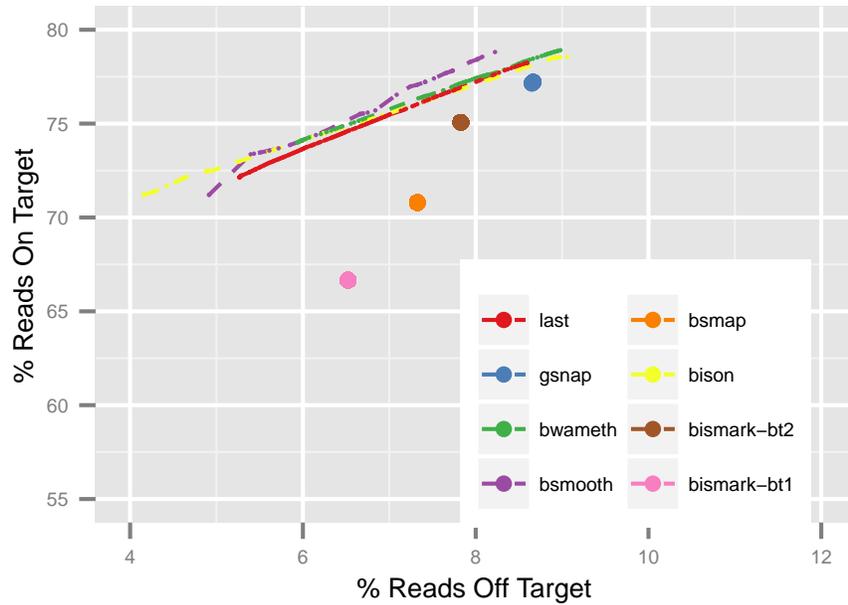


Figure 2: Percent of paired-end, *trimmed* 100-base reads on (y) and off (x) target for the tested aligners. Aligners that report mapping quality are shown as connected dots for each quality cut-off. Reads are limited to those considered as primary, mapped alignments by the aligner.

### 3.1 Resources

Table 2: Resources on real data

trimmed	program	CPU time(min)	mem(GB)	dataset
no	bis1	588.68	8.30	real
no	bis2	2067.23	9.88	real
no	bison	1639.25	10.34	real
no	bsmap	8498.52	23.06	real
no	bsmooth	5469.96	9.94	real
no	bwa	490.37	17.70	real
no	gsnap	5309.48	12.73	real
no	last	611.22	34.53	real
yes	bis1	581.60	8.28	real
yes	bis2	1921.00	9.64	real
yes	bison	1433.06	10.06	real
yes	bsmap	7259.29	23.19	real
yes	bsmooth	5189.16	9.46	real
yes	bwa	413.69	17.65	real
yes	gsnap	4565.57	13.93	real
yes	last	560.53	34.53	real

## 4 Comparison On Simulated Data

Paired-end, 100-base reads were simulated using the tool *Sherman* (<http://www.bioinformatics.babraham.a>) with 95% of the reads simulated from *mm10* and 5% of the reads from *E. coli*. Reads were aligned with the parameters in Supplemental Table 1. Reads from *E. coli* aligning to *mm10* were considered as off-target reads. The Y-axis in the plot is the percent of the *mm10* reads. Exact parameters sent to Sherman are here: <https://github.com/brentp/bwa-meth/blob/master/compare/src/gen-simulated.sh>

Note that these reads have no insertions or deletions. We simulated the *E. coli* reads to be 98% unmethylated and we used a per-base error rate of 1% for the simulations.

Supplemental Figure 3 below shows the result for these simulations.

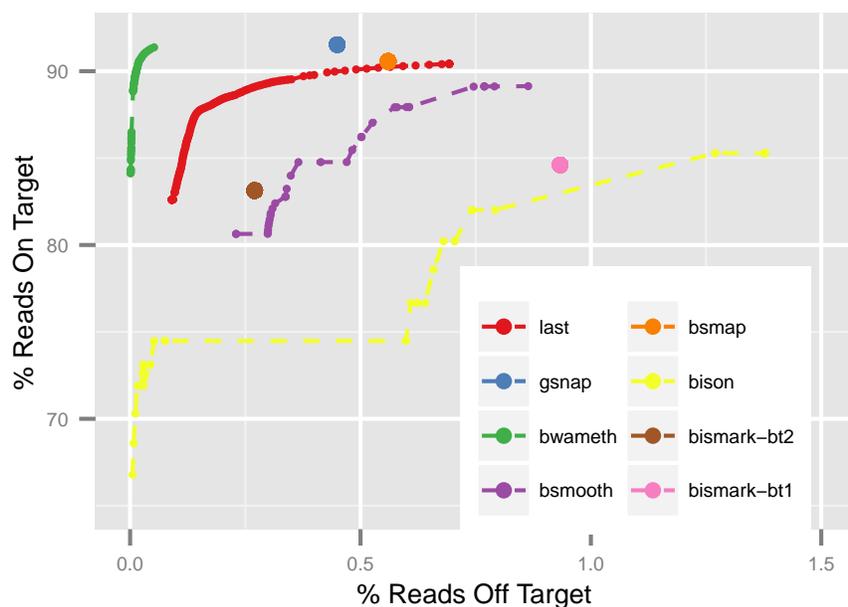


Figure 3: Percent of paired-end, *simulated* 100-base reads on (y) and off (x) target for the tested aligners. Aligners that report mapping quality are shown as connected dots for each quality cut-off. Reads are limited to those considered as primary, mapped alignments by the aligner.

Supplemental Figure 4 below shows the result for trimmed and simulated reads. Trimming removes some read-pairs if one or both reads had low-quality. We report the percent of reads aligned relative to the original, un-trimmed number from *mm10* since we are interested in the overall mapping rate.

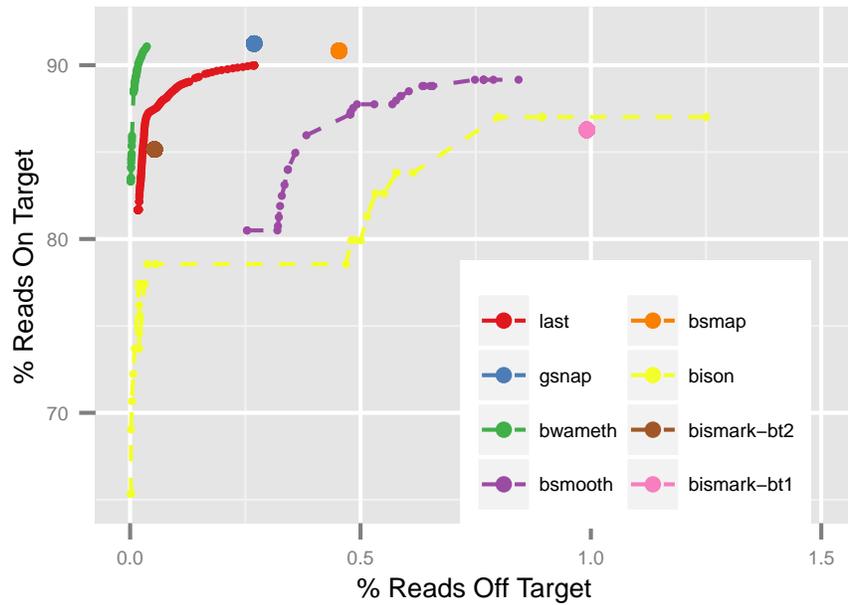


Figure 4: Percent of paired-end, *trimmed*, *simulated* 100-base reads on (y) and off (x) target for the tested aligners. Aligners that report mapping quality are shown as connected dots for each quality cut-off. Reads are limited to those considered as primary, mapped alignments by the aligner.

## 4.1 Resources

Table 3: Resources on simulated data with 1% error-rate (programs run with more processes will use more memory)

trimmed	program	CPU time(min)	mem(GB)	dataset
no	bis1	69.07	8.28	sim
no	bis2	271.13	9.66	sim
no	bison	222.23	9.51	sim
no	bsmap	227.35	22.99	sim
no	bsmooth	501.97	9.52	sim
no	bwa	145.56	19.40	sim
no	gsnap	265.41	10.84	sim
no	last	138.50	25.90	sim
yes	bis1	62.92	8.27	sim
yes	bis2	227.30	9.47	sim
yes	bison	154.81	9.16	sim
yes	bsmap	93.70	22.84	sim
yes	bsmooth	461.71	9.31	sim
yes	bwa	199.15	22.25	sim
yes	gsnap	91.98	10.78	sim
yes	last	113.80	25.90	sim

## 5 Comparison On Error-Free Simulated Data

Paired-end, error-free 100-base reads were simulated as above but without contamination or errors. Exact parameters sent to Sherman are here: <https://github.com/brentp/bwa-meth/blob/master/compare/src/gen-simulated.sh>

Supplemental Figure 5 below shows the result for these simulations.



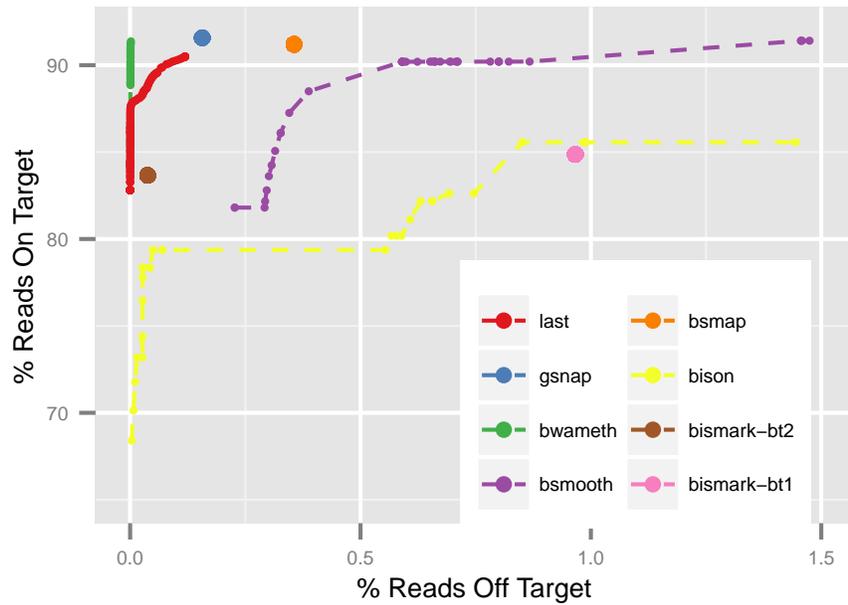


Figure 6: Percent of paired-end, *trimmed*, *error-free*, *simulated* 100-base reads on (y) and off (x) target for the tested aligners. Aligners that report mapping quality are shown as connected dots for each quality cut-off. Reads are limited to those considered as primary, mapped alignments by the aligner.

## 5.1 Resources

Table 4: Resources on error-free simulated data (programs run with more processes will use more memory)

trimmed	program	CPU time(min)	mem(GB)	dataset
no	bis1	61.05	8.27	error-free
no	bis2	289.67	9.46	error-free
no	bison	179.17	9.27	error-free
no	bsmap	66.03	22.72	error-free
no	bsmooth	492.62	9.37	error-free
no	bwa	170.66	21.11	error-free
no	gsnap	48.36	10.25	error-free
no	last	214.68	25.90	error-free
yes	bis1	64.32	8.27	error-free
yes	bis2	245.42	9.48	error-free
yes	bison	184.55	9.10	error-free
yes	bsmap	67.68	22.94	error-free
yes	bsmooth	481.99	9.38	error-free
yes	bwa	173.93	21.43	error-free
yes	gsnap	51.11	11.69	error-free
yes	last	138.03	25.90	error-free

## 6 Improved Accuracy for Stranded Capture Experiments

Below we show a reduction in the percent of off-target reads by considering only the strand targetted by our capture protocol.

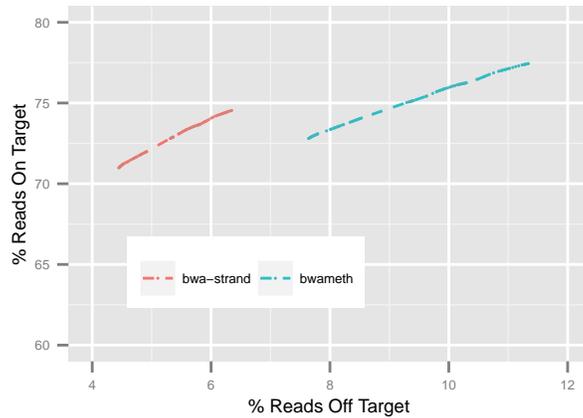


Figure 7: For a capture experiment that targets only one strand, we can improve accuracy by considering only reads that map to that strand. Here, we show the effect from the original `bwameth` alignments as shown in Figure 1 compared to the alignments restricted to the strand targeted by the capture method.

This is implemented in `bwa-meth` via the `-set-as-failed` flag so that reads mapping to a given strand are given the SAM flag indicating that they failed vendor quality control checks. Invocation would look like this:

```
bwameth.py \
  --reference $REF \
  --set-as-failed f \
  $READS_R1 $READS_R2
```

to set reads mapping to the forward, or original top strand as failed.

## 7 Mapping and Trimming

We note that trimming improves accuracy for most aligners but the difference is very small for `bwa-meth`.

## 8 `bwa-meth` Installation And Requirements

`Bwa-meth` depends on `samtools` and a single python library, `toolshed`. The latter can be installed by running `sudo python setup.py install` from the main directory of the `bwa-meth` project. `Samtools` is a C library installed on most systems and available at <https://github.com/samtools/samtools>. Further installation instructions are available in the `bwa-meth` README at <https://github.com/brentp/bwa-meth/>

For tabulation of methylation by CpG, `Bis-SNP` [?] is required. The java `.jar` file is available from: <http://sourceforge.net/projects/bissnp/files/>

For CNV detection from BS-Seq data, the R package *cn.mops* [?] is required. It is available from bioconductor [?] at: <http://bioconductor.org/packages/devel/bioc/html/cn.mops.html>

## 9 Additional Features

*Bwa-meth* and GSNAP are the only programs that output a BAM file that passes picard's ValidateSam without errors. Last does not report the proper pair information in all cases and none of the other aligners add a read-group. For the comparison, we added sorting and forced SAM output by the other aligners regardless of their default. *Bwa-meth* outputs a read-group for each sample by default and allows that to be customized.

*bwa-meth* can calculate bias of methylation estimates by location in the read:

```
python bias-plot.py input.bam ref.fa
```

This will create a bias-plot showing bases in the reads that should not be considered for methylation calls.

*Bwa-meth* defers tabulation of the methylation scores to Bis-SNP [?] by offering a simplified interface:

```
bwameth.py tabulate \  
  --trim 3,3 \  
  --map-q 30 \  
  --bissnp BisSNP-0.82.2.jar \  
  --reference /path/to/ref.fasta \  
  input.bam
```

Where the arguments are sent to Bis-SNP to, for example trim the first and last 3 bases from each read to avoid bias.

A full example on real data is at: <https://github.com/brentp/bwa-meth/tree/master/example/>