

A Poisson-Gaussian Denoising Dataset with Real Fluorescence Microscopy Images

Yide Zhang^{*}, Yinhao Zhu,^{*} Evan Nichols, Qingfei Wang, Siyuan Zhang, Cody Smith, Scott Howard
University of Notre Dame
Notre Dame, IN 46556, USA

{yzhang34, yzhu10, enichol3, qwang9, szhang8, csmith67, showard}@nd.edu

Abstract

Fluorescence microscopy has enabled a dramatic development in modern biology. Due to its inherently weak signal, fluorescence microscopy is not only much noisier than photography, but also presented with Poisson-Gaussian noise where Poisson noise, or shot noise, is the dominating noise source, instead of Gaussian noise that dominates in photography. To get clean fluorescence microscopy images, it is highly desirable to have effective denoising algorithms and datasets that are specifically designed to denoise fluorescence microscopy images. While such algorithms exist, there are no such datasets available. In this paper, we fill this gap by constructing a dataset - the Fluorescence Microscopy Denoising (FMD) dataset - that is dedicated to Poisson-Gaussian denoising. The dataset consists 12,000 real fluorescence microscopy images obtained with commercial confocal, two-photon, and wide-field microscopes and representative biological samples such as cells, zebrafish, and mouse brain tissues. We use imaging averaging to effectively obtain ground truth images and 60,000 noisy images with different noise levels. We use this dataset to benchmark 10 representative denoising algorithms and find that deep learning methods have the best performance. To our knowledge, this is the first microscopy image dataset for Poisson-Gaussian denoising purposes and it could be an important tool for high-quality, real-time denoising applications in biomedical research.

1. Introduction

Fluorescence microscopy is a powerful technique that permeates all of biomedical research [16]. Confocal [22], two-photon [10], and wide-field [25] microscopes are the most widely used fluorescence microscopy modalities that are vital to the development of modern biology. Fluorescence microscopy images, however, are inherently noisy

^{*}Equal contribution.

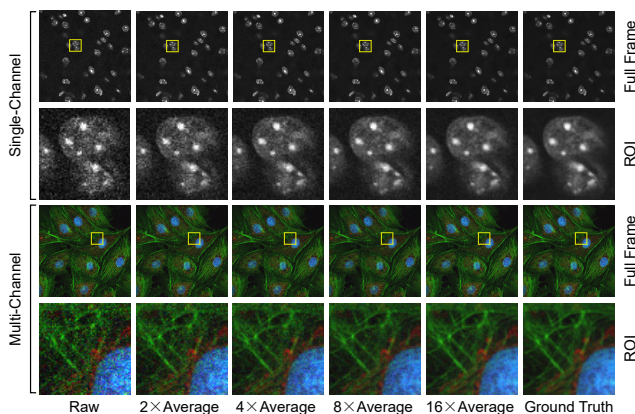


Figure 1. Examples of images with different noise levels and ground truth. The single-channel (gray) images are acquired with two-photon microscopy on fixed mouse brain tissues. The multi-channel (color) images are obtained with two-photon microscopy on fixed BPAE cells. The ground truth images are estimated by averaging 50 noisy raw images.

because the signal (photon flux) captured by a microscopic detector, such as a photomultiplier tube (PMT) or a charge coupled device (CCD) camera, is extremely weak compared to the signal in photography. And consequently, the measured optical signal in fluorescence microscopy is quantized due to the discrete nature of photons. One way to obtain clean images is to increase the power of the excitation laser or lamp, but the excitation power is not only limited by the dosages of light a biological sample can receive, but also fundamentally limited by the fluorescence saturation rate; i.e., the signal will stop to increase when the excitation power is too high [30]. Alternatively, one can get clean images by increasing the imaging time, e.g., pixel dwell time, exposure time, number of line or frame averages; this, however, may also cause photodamage to the sample. Moreover, for dynamic or real-time imaging, increasing the imaging time may be impossible since each image has to be captured within tens of milliseconds. Therefore, developing an

algorithm to effectively denoise (reduce the noisiness of) a fluorescence microscopy image is of great importance to biomedical research. Meanwhile, a high-quality denoising dataset is necessary to benchmark and evaluate the effectiveness of the denoising algorithm.

Most of the image denoising algorithms and datasets are created for Gaussian noise dominated images, with a recent focus on denoising with real noisy images, such as smart phones [1] or DSLR images [23]. However, there is a lack of a reliable Poisson noise dominated denoising dataset comprising of *real* fluorescence microscopy images. The goal of this work is to fill this gap. More specially, we create a Poisson-Gaussian denoising dataset - the Fluorescence Microscopy Denoising (FMD) dataset - consisting of 12,000 real noisy microscopy images which cover the three most widely used imaging modalities, i.e., confocal, two-photon, and wide-field, as well as three representative biological samples including cells, zebrafish, and mouse brain tissues. With high-quality commercial microscopy, we use image averaging to effectively obtain ground truth images and noisy images with five different noise levels. Some image averaging examples are shown in Figure 1. We further use this dataset to benchmark classic denoising algorithms and recent deep learning models, with or without ground truth. Our FMD dataset is publicly available¹, including a convenient PyTorch data loader. To our knowledge, this is the first dataset constructed from real noisy fluorescence microscopy images and designed for Poisson-Gaussian denoising purposes.

2. Related Work

There are consistent efforts in constructing denoising dataset with real images to better capture the real-world noise characteristics and evaluate denoising algorithms, such as RENOIR [4], Darmstadt Noise Dataset [23], Smartphone Image Denoising Dataset [1], and PolyU Dataset [27]. Those datasets contain real images taken from either DSLR or smartphones with different ISOs and different number of scenes. The dominating noise in those images is Gaussian or Poisson-Gaussian in real low-light conditions. However, there is no dedicated dataset for Poisson noise dominated images, which are inherently different from Gaussian denoising datasets. This work is dedicated for fluorescence microscopy denoising where the images are corrupted by Poisson-Gaussian noise; in particular, Poisson noise, or shot noise, is the dominant noise source.

Image averaging is the most used method to obtain ground truth images when constructing denoising dataset. The main efforts are spent on image pre-processing, such as image registration to remove the spatial misalignment of an

image sequence with the same FOV [3, 1], intensity scaling due to the changes of light strength or analog gain [23], and methods to cope with clipped pixels due to over exposure or low-light conditions [4]. The images captured by commercial microscopes in our dataset turns out to be well aligned, and the analog gain is carefully chosen to avoid clipping and to utilize the full dynamic range.

There are two main approaches to denoise an image corrupted by Poisson-Gaussian noise. One way is to directly apply an effective denoising algorithm, such as the PURE-LET method [18], which is designed to handle the Poisson-Gaussian denoising problem based on the statistics of the noise model. Another approach is using a nonlinear variance-stabilizing transformation (VST) to convert the Poisson-Gaussian denoising problem into a Gaussian noise removal problem, which is well studied with a considerable amount of effective denoising algorithms to choose from, such as NLM, BM3D, KSVD, EPLL, and WNNM [7, 9, 2, 31, 12] etc. The VST-based denoising process generally involves three steps. First, the noisy raw images are transformed using a VST designed for the noise model. In our case, we use the generalized Anscombe transformation (GAT) that is designed for Poisson-Gaussian noise [20]. The VST is able to remove the signal-dependency of the Poisson component, whose noise variance varies with the expected pixel value, and results in a modified image with signal-independent Gaussian noise only and a constant, i.e., unitary, noise variance. Next, a Gaussian denoising algorithm is applied to the transformed image. And finally, the Gaussian-denoised data is transformed back via an inverse VST algorithm, such as the exact unbiased inverse transformation [20], and the estimation of the noise-free image is obtained.

Recently there is an increasing interest in deep learning based methods for image denoising, where fully convolutional networks (FCNs) [17] are used for this image-to-image regression problem. With residual learning and batch normalization, DnCNN [28] reports better performance than traditional denoising methods such as BM3D. Further development towards blind image denoising includes incorporating non-uniform noise level map in the input of FFDNet [29], or noise estimation network as in CBDNet [13], or utilizing the non-local self-similarity in UDNNet [14] and [24]. These methods all require clean images to supervise the training. There are also progress on denoising methods without paired clean images [8] using generative adversarial networks to learn the noise model. In [15], a Noise2Noise model is trained without clean images at all and outperforms VST+BM3D by almost 2dB on synthetic Poisson noise.

We perform intensive study of the noise statistics of the FMD dataset and show the noise is indeed Poisson-dominated for two-photon and confocal microscopy, and

¹https://drive.google.com/open?id=1aygMzSDdoq63IqSk-ly8cMq0_owup8UM

has larger Gaussian component for wide-field microscopy. We then benchmark 10 representative denoising algorithms on the FMD dataset, and show better denoising performance with deep learning models than traditional methods on the real noisy images.

3. Noise Modeling in Fluorescence Microscopy

The microscopy imaging system is modeled with a Poisson-Gaussian noise model [11, 20]. The model is composed of a Poisson noise component that accounts for the signal-dependent uncertainty, i.e., shot noise, and an additive Gaussian noise component which represents the signal-independent uncertainty such as thermal noise. Specifically, let z_i , $i = 1, 2, \dots, N$, be the measured pixel values obtained with a PMT or a CCD, and

$$z_i = y_i + n_i = y_i + n_p(y_i) + n_g, \quad (1)$$

where y_i is the ground truth and n_i is the noise of the pixel; the noise n_i is composed of two mutually independent parts, n_p and n_g , where n_p is a signal-dependent Poisson noise component that is a function of y_i , and n_g is a signal-independent zero-mean Gaussian component. Denoting $a > 0$ as the conversion or scaling coefficient of the detector, i.e., a single detected photon corresponds to a measured pixel value of a , and $b \geq 0$ as the variance of the Gaussian noise, we can describe the Poisson and Gaussian (normal) distributions as

$$(y_i + n_p(y_i))/a \sim \mathcal{P}(y_i/a), \quad n_g \sim \mathcal{N}(0, b). \quad (2)$$

Note that a is related to the quantum efficiency of the detector. Assuming that the Poisson and Gaussian processes are independent, the probability distribution of z_i is the convolution of their individual distributions, i.e.,

$$p(z_i) = \sum_{k=0}^{+\infty} \left(\frac{\left(\frac{y_i}{a}\right)^k}{k!} e^{-\frac{y_i}{a}} \times \frac{1}{\sqrt{2\pi b}} e^{-\frac{(z_i - ak)^2}{2b}} \right). \quad (3)$$

The denoising problem of a microscopy image is then to estimate the underlying ground truth y_i given the noisy measurement of z_i .

To denoise a fluorescence microscopy image, one can use algorithms that are specifically designed for Poisson-Gaussian denoising. A more common approach is using VST to stabilize the variance such that the denoising task can be tackled by a well-studied Gaussian denoising method. As a representative VST method, GAT transforms the measured pixel value z_i in the image to

$$f(z_i) = \frac{2}{a} \sqrt{\max\left(az_i + \frac{3}{8}a^2 + b, 0\right)}, \quad (4)$$

which stabilizes its noise variance to approximately unity, i.e., $\text{Var}\{f(z_i)\} \approx 1$. A Gaussian denoising algorithm,

such as NLM and BM3D, can then be applied to $f(z_i)$ because its noise can be considered as a signal-independent Gaussian process with zero mean and unity variance. Once the denoised version of $f(z_i)$, $D(z_i)$ is obtained, an inverse VST is used to estimate the signal of interest y_i . However, simply applying an algebraic inverse f^{-1} to D will generally result in a biased estimate of y_i . An asymptotically unbiased inverse can mitigate the bias, but the denoising accuracy will be problematic for images with low signal levels, a common property of fluorescence microscopy images [6]. To address this problem, we use the exact unbiased inverse transformation, which can estimate the signal of interest accurately even at low signal levels [20]. Practically, since the exact unbiased inverse requires tabulation of parameters, one can employ a closed-form approximation of it [19], i.e.,

$$\tilde{I}(D) = \frac{1}{4}D^2 + \frac{1}{4}\sqrt{\frac{3}{2}}D^{-1} - \frac{11}{8}D^{-2} + \frac{5}{8}\sqrt{\frac{3}{2}}D^{-3} - \frac{1}{8}. \quad (5)$$

The closed-form approximation ensures the denoising accuracy while reducing the computational cost, and the estimated noise-free signal is $\tilde{y}_i = \tilde{I}[D(z_i)]$.

To evaluate and benchmark the performances of different denoising algorithms, a ground truth and images with various noise levels are needed, which can be obtained by averaging a series of noisy raw fluorescence microscopy images on the same FOV. The averaging is performed after ensuring that no image shift larger than a half-pixel can be detected by an image registration algorithm. Since for different raw images, their Poisson-Gaussian random processes are independent, the average of S noisy raw images, v_i^S , can be written as

$$\begin{aligned} v_i^S &= \frac{1}{S} \sum_{j=1}^S z_i^j = \frac{a}{S} \sum_{j=1}^S \frac{y_i + n_p^j(y_i)}{a} + \frac{1}{S} \sum_{j=1}^S n_g^j \\ &\sim \frac{a}{S} \mathcal{P}\left(\frac{Sy_i}{a}\right) + \frac{1}{S} \mathcal{N}(0, Sb), \end{aligned} \quad (6)$$

where n_p^j and n_g^j are the noise realizations of the j -th noisy image. Based on the properties of Poisson and Gaussian distributions, the mean and variance of the averaged image, v_i^S , can be written as

$$\mathbb{E}[v_i^S] = y_i, \quad \text{Var}[v_i^S] = \frac{a}{S}y_i + \frac{b}{S}. \quad (7)$$

As the number of noisy images used for averaging increases, the noise of ground truth estimation, $\sqrt{\text{Var}[v_i^S]}$, decreases, while the ground truth signal, $\mathbb{E}[v_i^S]$ is invariant; therefore, image averaging is equivalent to increasing the signal-to-noise ratio (SNR) of estimating the ground truth. We make $S = 1, 2, 4, 8, 16$ to create images with five different noise levels, and $S = 50$ to generate the ground truth.

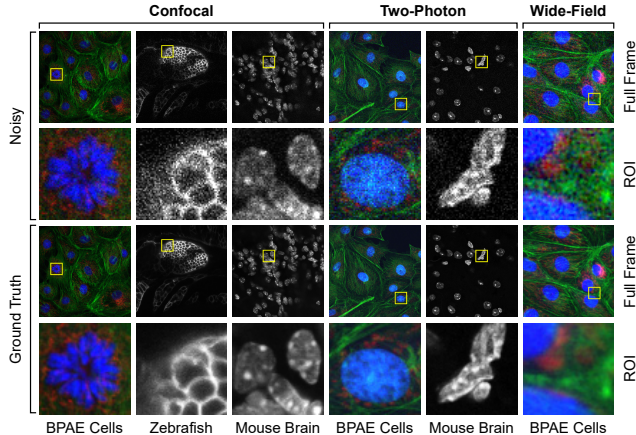


Figure 2. Examples of raw fluorescence microscopy images and their estimated ground truth from our FMD dataset. Shown here are FOVs from different microscopy modalities on different biological samples.

As demonstrated in [3] and also shown in Sec. 4.3, for fluorescence microscopy images, little image quality improvement can be seen after including around 40 images in averaging.

4. Dataset

In this section, we describe the experimental setup that we used to acquire fluorescence microscopy images. We then discuss how the raw images are utilized to estimate ground truth as well as images with different noise levels. Finally we present the statistics as well as the estimated noise levels of our dataset.

4.1. Image Acquisition Setup

Our FMD dataset covers the three main modalities of fluorescence microscopy: confocal, two-photon, and wide-field. All images were acquired with high-quality commercial fluorescence microscopes and imaged with real biological samples, including fixed bovine pulmonary artery endothelial (BPAE) cells [labeled with MitoTracker Red CMXRos (mitochondria), Alexa Fluor 488 phalloidin (F-actin), and DAPI (nuclei); Invitrogen FluoCells F36924], fixed mouse brain tissues (stained with DAPI and cleared), and fixed zebrafish embryos [EGFP labeled *Tg(sox10:megfp)* zebrafish at 2 days post fertilization]. All animal studies were approved by the university’s Institutional Animal Care and Use Committee. To acquire noisy microscopy images for denoising purposes, we kept an excitation laser/lamp power as low as possible for all imaging modalities. Specifically, the excitation power was low enough to generate a very noisy image, yet high enough such that the image features were discernible. We also manually set the detector/camera gain to a proper value to avoid

clipping and to fully utilize the dynamic range.

The confocal and two-photon images were acquired with a Nikon A1R-MP laser scanning confocal microscope equipped with a Nikon Apo LWD 40 \times , 1.15 NA water-immersion objective. The confocal and two-photon images were 512 \times 512 pixels with a pixel size of 300 nm and a pixel dwell time of 2 μ s. The A1R-MP microscope has multiple detectors (PMTs) in parallel, so for multi-channel (color) fluorescence imaging with the BPAE cells, all three images were acquired simultaneously. For confocal imaging, the excitation was generated by a LU4/LU4A laser unit, the pinhole size was set to 1.2 Airy unit, and the imaging conditions for different samples were as follows: BPAE nuclei, 405 nm excitation, 0.5% laser power, 110 PMT gain; BPAE F-actin, 488 nm excitation, 0.5% laser power, 110 PMT gain; BPAE mitochondria, 561 nm excitation, 0.5% laser power, 110 PMT gain; mouse brain, 405 nm excitation, 0.5% laser power, 115 PMT gain; zebrafish embryo, 488 nm excitation, 10% laser power, 140 PMT gain. For two-photon microscopy, the excitation was generated by a Spectra-Physics Mai Tai DeepSee femtosecond laser, and for all two-photon images, the laser power was set to 0.5%, the PMT gain to 130, and the excitation wavelength to 780 nm. Note that our dataset did not include two-photon images of the zebrafish sample because during two-photon imaging, very strong two-photon auto-fluorescence signals from the zebrafish were observed, which severely degraded the imaging quality.

The wide-field images were acquired with a Nikon Eclipse 90i wide-field fluorescence microscope equipped with a Nikon Plan Fluor 40 \times , 0.75 NA objective. The excitation was generated by a halogen lamp (with ND16 neutral-density filter) and the images were captured by a DS-Fi1-U2 camera with an exposure time of 200 ms and a gain of 46. The raw image size was 1280 \times 960 and the pixel size was 170 nm. These images were cropped to 512 \times 512 before being processed for our dataset. Note that our dataset only covered wide-field images of the BPAE cells because wide-field microscopy could not image well in animal tissues such as mouse brain and zebrafish embryo, where strong out-of-focus fluorescence would blur out the wide-field images. Since the BPAE cells were stained with three different fluorophores while only one detector (CCD camera) was available in the 90i microscope, we imaged three times for the same FOV, each time with a different filter block (DAPI for nuclei, FITC for F-actin, TRITC for mitochondria), to acquire the multi-channel (color) fluorescence image of the cells.

For any imaging modality, each sample was imaged with 20 different FOVs, and each FOV was repeatedly captured for 50 times as 50 noise realizations. The acquired images were preprocessed and used for noisy image and ground truth estimation as described in Sec. 4.2. Figure 2 shows

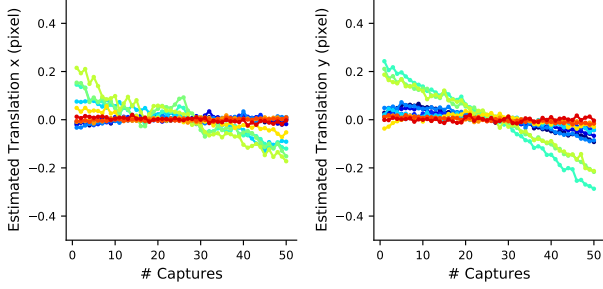


Figure 3. Estimated translation along x and y axes, both within a half-pixel (0.5). The estimation is performed on the 20-th FOV of each imaging configuration. Each line in a plot shows the estimation of one of the 12 configurations (different modalities on different samples).

some example images of a single FOV from different imaging modalities and different samples.

4.2. Noisy Image and Ground Truth Estimation

Image registration The approach to estimate ground truth by averaging a sequence of captures usually comes with the issue of spatial misalignment, which is typical in photos taken by smartphones and DSLR. We use intensity-based image registration to register a sequence of image with the same FOV against the mean image of the sequence, but find that the estimated global translations in both x and y axis are less than a half-pixel (0.5), as shown in Figure 3. Translation in sub-pixel smooths out noisy images, thus destroys the realness of Poisson noise which is the main characteristic of our dataset. In short, the image sequence obtained by the commercial fluorescence microscope is already well aligned, thus image registration is not performed.

Different noise levels As described in Section 4.1, the raw images are acquired with a low excitation power thus a relatively high noise level (low SNR) to increase the difficulty of denoising task. Meanwhile, the raw images with high noise levels allows us to create images with lower noise levels by image averaging. Particularly, we obtain averaged images with four extra noise levels by averaging S ($S = 2, 4, 8, 16$) raw images, respectively, within the same sequence (FOV) of 50. We sequentially select each image within the sequence; for each selected image, $S - 1$ images next to it are circularly selected; the S selected images in total are used for averaging. Using this circular averaging method, we are able to obtain the same number of averaged images as the number of raw images in the sequence, i.e., 50; meanwhile, the newly generated 50 raw images can be considered as 50 different noise realizations. In this way, the amount of noisy images in the dataset can be increased to five-fold ($S = 1, 2, 4, 8, 16$). Some example images with different noise levels are shown in Figure 1. As also shown

Modality	Samples	a	b
CF	BPAE (Nuclei)	1.39×10^{-2}	-2.16×10^{-4}
CF	BPAE (F-actin)	1.37×10^{-2}	-1.85×10^{-4}
CF	BPAE (Mito)	1.21×10^{-2}	-1.54×10^{-4}
CF	Zebrafish	9.43×10^{-2}	-1.60×10^{-3}
CF	Mouse Brain	1.94×10^{-2}	-2.68×10^{-4}
TP	BPAE (Nuclei)	3.31×10^{-2}	-8.39×10^{-4}
TP	BPAE (F-actin)	2.55×10^{-2}	-5.43×10^{-4}
TP	BPAE (Mito)	2.10×10^{-2}	-4.57×10^{-4}
TP	Mouse Brain	3.38×10^{-2}	-9.16×10^{-4}
WF	BPAE (Nuclei)	2.29×10^{-4}	2.35×10^{-4}
WF	BPAE (F-actin)	1.94×10^{-3}	1.91×10^{-4}
WF	BPAE (Mito)	3.55×10^{-4}	1.95×10^{-4}

Table 1. Estimation of noise parameters (a , b) of the FMD dataset. The shown a and b are average estimation values of 20 raw noisy images from 20 different FOVs (a raw image from each FOV). CF, confocal; TP: two-photon; WF: wide-field.

in Table 2, the peak signal-to-noise ratio (PSNR) of the averaged images increases as the number of raw images used for averaging increases.

Ground truth estimation We estimate the ground truth by averaging all 50 captures of the same FOV, similar to the approaches employed in [3] and [18]; hence in the FMD dataset, each FOV has only one ground truth that is shared by all noise realizations from this FOV. As demonstrated in [3] and also shown in Sec. 4.3, the image quality or noise characteristics of a fluorescence microscopy image will see little improvement after including around 40 images in the average; therefore we choose 50 captures as our criterion to obtain the ground truth. As shown in Equations (6) and (7), the ground truth y_i for images with different noise levels z_i^j is the same, and image averaging is equivalent to sampling from a Poisson-Gaussian distribution with a higher SNR. Regardless of the number of images used for averaging, the mean stays the same and equals to the ground truth. Figure 1 shows two ground truth images as well as their corresponding noise realizations.

4.3. Dataset Statistics and Noise Estimation

Taking the combination of each sample (the BPAE cells are considered as three samples due to its fluorophore composition) and each microscopy modality as a configuration, the FMD dataset includes 12 different imaging configurations that are representative of almost all fluorescence microscopy applications in practice. For each configuration, we capture 20 different FOVs of the sample, and for each FOV, we acquire 50 raw images. Meanwhile, the 50 raw images in a FOV can be extended to five-fold using the circular averaging method described in Sec. 4.2. Therefore, in total, the dataset has $12 \times 20 = 240$ FOVs or ground truth images, $240 \times 50 = 12,000$ raw images, and $12,000 \times 5 = 60,000$

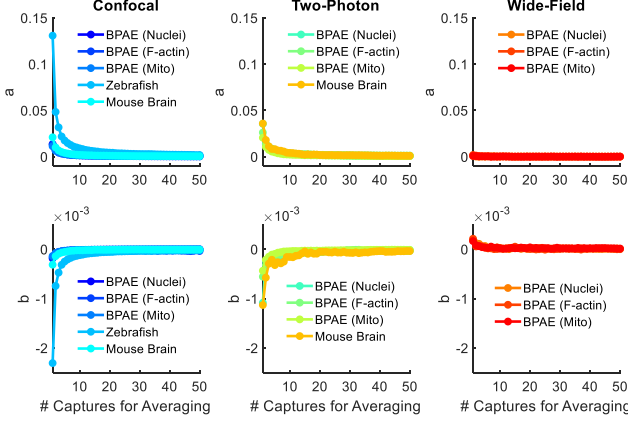


Figure 4. Estimated noise parameters (a and b) of averaged images obtained with different raw image numbers in the average. The estimation is performed on the second FOV of each imaging configuration.

noisy images as noise realizations.

While there are blind denoising methods (e.g., DnCNN) that are able to denoise an image without any additional information, most denoising algorithms such as NLM and BM3D, however, require an estimate of the noise levels presented in the image. In this work, we employ the noise estimation method in [11] to estimate the Poisson-Gaussian noise parameters, a and b , described in Sec. 3. The estimated values of a and b not only are needed in the benchmark of various denoising algorithms, they also reflect the characteristics of the noise presented in an images. Specifically, since Poisson-Gaussian noise is a mixture of both Poisson and Gaussian noises, which are parameterized by a and b , respectively, an image with a large estimate value of a but a small b may be considered as a Poisson noise dominated image, while a small a with a large b can indicate that the image is Gaussian noise dominated. In fluorescence microscopy, however, it is unlikely to have a Gaussian noise dominated image due to the low signal levels; most fluorescence microscopy images are Poisson noise, or shot noise, dominated, with certain types of microscopes, such as wide-field ones, have a considerable amount of Gaussian noise involved [5, 21]. Note that the noise estimation program from [11] could generate a negative b value when the Gaussian noise component is very small compared to the Poisson component. This, however, does not mean that the image has a “negative” Gaussian noise variance. More details can be found in [11].

We evaluate the noise characteristics of our FMD dataset by estimating the noise parameters of raw noisy image (1 in each FOV, 240 in total). The estimated a and b are then grouped according to their corresponding imaging configurations (20 FOVs in each configuration, 12 configurations in total) and averaged. The results are presented in Ta-

ble 1. For confocal and two-photon microscopy, the estimated a are comparably large while the b are negative; hence confocal and two-photon images are Poisson noise dominated. For wide-field microscopy, however, the a are much smaller than above, possibly due to the much lower sensitivity of CCD cameras used in wide-field microscopy compared to the PMTs used in confocal and two-photon microscopy; meanwhile, the b are now all positive, which indicates that wide-field images have a mixed Poisson-Gaussian noise with a considerable amount of Gaussian noise presented. We further evaluate the effect of image averaging on its noise characteristics. Figure 4 shows the estimated a and b values when different number of images, S , are included in the average. The results are in good agreement with the theory in Eq. 7 and the observations in Table 1, as the estimated parameters follows the trend of a/S and b/S , and their initial values ($S = 1$) are close to the ones in Table 1. Figure 4 also shows that the values of a and b exhibit little change when the number of captures used for averaging is more than 40; this confirms the observation reported in [3] that the image quality or noise characteristics of a fluorescence microscopy image will see little improvement after including around 40 images in the average.

5. Benchmark

In this section we benchmark several representative denoising methods, including deep learning models, on our fluorescence microscopy images with real Poisson-Gaussian noise. We show that deep learning models perform better than traditional methods on the FMD dataset.

5.1. Setup

The FMD dataset is split to training and test sets, where the test set is composed of images randomly selected from the 19-th FOV of each imaging configuration (the rest 19 FOVs are for training and validation purposes). The mixed test set consists of 4 images randomly selected from the 19-th FOV of 12 imaging configurations (combination of microscopy modalities and biological samples), organized in different noise levels. Thus we have 5 mixed test sets each of which have 48 noisy images with a specific noise level corresponding to 1 (raw), 2, 4, 8, or 16 times averaging. We also test the denoising algorithms on all 50 images from the same FOV of a specific imaging configuration, also organized in different noise levels, with denoising results shown in supplementary materials.

Considering GPU memory constraint for training fully convolutional networks [28, 15] on large images, we crop the raw images of size 512×512 to four non-overlapping patches of size 256×256 . We evaluate the computation time on single core of Intel Xeon CPU E5-2680, and additionally on single core of Nvidia GeForce GTX 1080 Ti GPU for deep learning models.

The 10 benchmarked algorithms in this work can be divided into three categories. The first category is for the methods that are specifically designed for Poisson-Gaussian denoising; we benchmark PURE-LET [18], an effective and representative Poisson-Gaussian denoising algorithm. The second category is for using well-studied Gaussian denoising methods in combination with VST and inverse VST; we combine GAT and the exact unbiased inverse transformation with classical denoising algorithms including NLM [7], BM3D [9], KSVD and its two variants KSVD(D) (overcomplete DCT dictionary) and KSVD(G) (global or given dictionary) [2], EPLL [31], and WNNM [12]. The last category is for deep learning based methods; we benchmark DnCNN [28] and Noise2Noise [15]. Note that the estimation of noise parameters a and b are required for the algorithms in the first and second categories to work. The estimation is performed according to Sec. 4.3 and then the images as well as the estimated parameters are sent to the denoising algorithms.

For benchmarking deep learning methods, unlike previous work [1] that directly tests with the pre-trained models, we re-train these models with the same network architecture and similar hyper-parameters on the FMD dataset from scratch. Specifically, we compare two representative models, one of which requires ground truth (DnCNN) and the other does not (Noise2Noise).

5.2. Results and Discussion

The benchmark denoising results on the mixed test set is shown in Table 2, including PSNR, structural similarity index (SSIM) [26] and denoising time. From the table, BM3D (in combination with VST) is still the most versatile traditional denoising algorithm regarding its high PSNR and relatively fast denoising speed. PURE-LET, though its PSNR is not the highest, is the fastest denoising method among all the benchmarked algorithms thanks to its specific design for Poisson-Gaussian denoising. Finally, deep learning models outperform the other 8 methods by a significant margin in all noise levels, both in terms of PSNR and SSIM, even though they are *blind* to noise levels. This is different from the observation made before in [1, 23], probably because the nature of Poisson dominated noise is different from Gaussian noise while most of the denoising methods are developed for Gaussian noise model. Even if we applied the VST before Gaussian denoising, the transformed noise may still be different from a pure Gaussian one. More importantly, here the models are re-trained with our FMD dataset instead of pre-trained on other datasets.

For the Noise2Noise model, we add a batch normalization layer after every convolution layer as well as a Tanh activation function in the output layer. We find that this modification helps stabilize the training process. The training data includes all imaging configurations and noise lev-

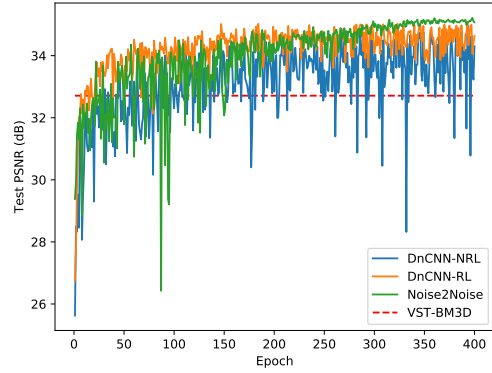


Figure 5. Test PSNR on the mixed test set with raw images during training. Each training epoch contains 4560 ($= 5 \times 12 \times 19 \times 4$) images of size 256×256 . Given enough training time (e.g. 400 epochs), Noise2Noise eventually outperforms DnCNN and VST-BM3D. Also notice for DnCNN, residual learning does help improve denoising.

els; thus we use one trained model to perform blind denoising on various imaging configurations and noises. We confirm that overall the Noise2Noise model has similar denoising performance as DnCNN, but without the need of clean images, and with almost 2dB higher than VST+BM3D in PSNR [15]. It even performs slightly better than DnCNN in the high noise domain, which is desirable in practice. We also train DnCNN without residual learning (DnCNN-NRL) where the model directly outputs the denoised image instead of the residual between clean and noisy images, and confirm it is worse than the model with residual learning (DnCNN-RL), as has been reported in [28]. The test performance for the mixed test set with raw images during training is shown in Figure 5. We also show benchmark results of the 10 algorithms on single-channel (gray) and multi-channel (color) confocal images in Figures 6 and 7, respectively, where the PSNR and SSIM of the color images are the mean values of that of their three channels.

Regarding the denoising speed, deep learning models have similar denoising time with that of VST-BM3D, and are even slower than PURE-LET, when running on CPU. However, the denoising time can be reduced to less than 1 ms when running on GPU, which potentially enables real-time denoising up to 100 frames per second, which is out of reach of traditional denoising methods. With such a denoising speed and high performance, deep learning denoising methods could dramatically benefit real-time fluorescence microscopy imaging, which allows biomedical researchers to observe the fast and dynamic biological processes in a much improved quality and to see processes that cannot be clearly seen before.

Methods	Number of raw images for averaging					Time
	1	2	4	8	16	
Raw	27.22 / 0.5442	30.08 / 0.6800	32.86 / 0.7981	36.03 / 0.8892	39.70 / 0.9487	-
VST+NLM [7]	31.25 / 0.7503	32.85 / 0.8116	34.92 / 0.8763	37.09 / 0.9208	40.04 / 0.9540	137.10 s
VST+BM3D [20]	32.71 / 0.7922	34.09 / 0.8430	36.05 / 0.8970	38.01 / 0.9336	40.61 / 0.9598	5.67 s
VST+KSVD [2]	32.02 / 0.7746	33.69 / 0.8327	35.84 / 0.8933	37.79 / 0.9314	40.36 / 0.9585	341.21 s
VST+KSVD(D) [2]	31.77 / 0.7712	33.45 / 0.8292	35.67 / 0.8908	37.69 / 0.9300	40.32 / 0.9579	67.96 s
VST+KSVD(G) [2]	31.98 / 0.7752	33.64 / 0.8327	35.83 / 0.8930	37.82 / 0.9312	40.44 / 0.9584	58.82 s
VST+EPLL [31]	32.61 / 0.7876	34.07 / 0.8414	36.08 / 0.8970	38.12 / 0.9349	40.83 / 0.9618	288.63 s
VST+WNNM [12]	32.52 / 0.7880	34.04 / 0.8419	36.04 / 0.8973	37.95 / 0.9334	40.45 / 0.9587	451.89 s
PURE-LET [18]	31.95 / 0.7664	33.49 / 0.8270	35.29 / 0.8814	37.25 / 0.9212	39.59 / 0.9450	2.61 s
DnCNN [28]	34.73 / 0.9043	36.13 / 0.9281	37.64 / 0.9476	39.34 / 0.9597	41.72 / 0.9725	35.05 s
Noise2Noise [15]	35.19 / 0.9117	36.16 / 0.9284	37.68 / 0.9471	39.26 / 0.9585	41.28 / 0.9705	9.43 s

Table 2. Denoising performance using the mixed test set, which includes confocal, two-photon, and wide-field microscopy images. PSNR (dB), SSIM, and denoising time (seconds) are obtained by averaging over 48 noise realizations in the mixed test set for each of 5 noise levels. Results of DnCNN and Noise2Noise are obtained by training on dataset with all noise levels. All 50 captures of each FOV (except the 19-th FOV which is reversed for test) are included in the training set, with 1 (DnCNN) or 2 (Noise2Noise) samples of which randomly selected from each FOV when forming mini-batches during training for 400 epochs. Note that test time for deep learning models on GPU is faster in orders of magnitude, i.e. 0.62 ms for DnCNN and 0.99 ms for Noise2Noise on single GPU in our experiment.

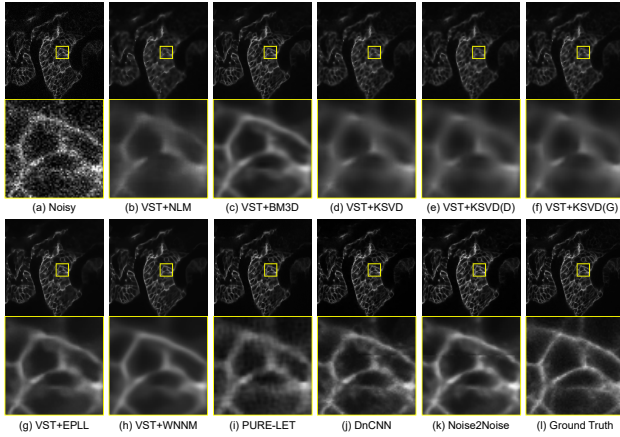


Figure 6. Benchmark results for single-channel (gray) images (zebrafish embryo under confocal microscopy). (a) Noisy, PSNR: 22.71 dB, SSIM: 0.4441. (b) VST+NLM, PSNR: 28.49 dB, SSIM: 0.7952. (c) VST+BM3D, PSNR: 31.99 dB, SSIM: 0.8862. (d) VST+KSVD, PSNR: 29.25 dB, SSIM: 0.8234. (e) VST+KSVD(D), PSNR: 29.04 dB, SSIM: 0.8212. (f) VST+KSVD(G), PSNR: 29.23 dB, SSIM: 0.8232. (g) VST+EPLL, PSNR: 31.71 dB, SSIM: 0.8711. (h) VST+WNNM, PSNR: 31.22 dB, SSIM: 0.8702. (i) PURE-LET, PSNR: 30.59 dB, SSIM: 0.8332. (j) DnCNN, PSNR: 32.27 dB, SSIM: 0.9010. (k) Noise2Noise, **PSNR: 32.85 dB, SSIM: 0.9070**. (l) Ground Truth.

6. Conclusion

In this work, we have constructed a dedicated denoising dataset of real fluorescence microscopy images with Poisson-Gaussian noise, which covers most microscopy modalities. We have used image averaging to obtain ground truth and noisy images with 5 different noise levels. With

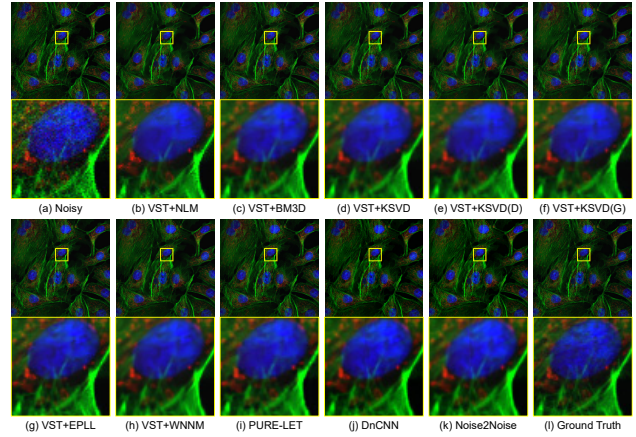


Figure 7. Benchmark results for multi-channel (color) images (BPAE cells under confocal microscopy). (a) Noisy, PSNR: 30.67 dB, SSIM: 0.7902. (b) VST+NLM, PSNR: 34.74 dB, SSIM: 0.9108. (c) VST+BM3D, PSNR: 35.86 dB, SSIM: 0.9338. (d) VST+KSVD, PSNR: 35.72 dB, SSIM: 0.9209. (e) VST+KSVD(D), PSNR: 35.47 dB, SSIM: 0.9139. (f) VST+KSVD(G), PSNR: 35.63 dB, SSIM: 0.9176. (g) VST+EPLL, PSNR: 35.72 dB, SSIM: 0.9335. (h) VST+WNNM, PSNR: 35.89 dB, SSIM: 0.9322. (i) PURE-LET, PSNR: 35.18 dB, SSIM: 0.9262. (j) DnCNN, PSNR: 35.95 dB, **SSIM: 0.9432**. (k) Noise2Noise, **PSNR: 36.45 dB, SSIM: 0.9428**. (l) Ground Truth.

this dataset, we have benchmarked representative denoising algorithms for Poisson-Gaussian noise including the most recent deep learning models. The benchmark results show that deep learning denoising models trained on our FMD dataset outperforms other methods by a large margin across all imaging modalities and noise levels. We have made our FMD dataset publicly available as a benchmark for Poisson-

Gaussian denoising research, which, we believe, will be especially useful for researchers that are interested in improving the imaging quality of fluorescence microscopy.

Acknowledgments

This material is based upon work supported by the National Science Foundation under Grant No. CBET-1554516. Yide Zhang's research was supported by the Berry Family Foundation Graduate Fellowship of Advanced Diagnostics & Therapeutics (AD&T), University of Notre Dame. The authors further acknowledge the Notre Dame Integrated Imaging Facility (NDIIF) for the use of the Nikon A1R-MP confocal microscope and Nikon Eclipse 90i wide-field microscope in NDIIF's Optical Microscopy Core.

References

- [1] A. Abdelhamed, S. Lin, and M. S. Brown. A high-quality denoising dataset for smartphone cameras. In *Proceedings of the IEEE Conference on Computer Vision and Pattern Recognition*, pages 1692–1700, 2018.
- [2] M. Aharon, M. Elad, and A. Bruckstein. K-SVD: An Algorithm for Designing Overcomplete Dictionaries for Sparse Representation. *IEEE Transactions on Signal Processing*, 54(11):4311–4322, 2006.
- [3] N. S. Alexander, G. Palczewska, P. Stremplewski, M. Wojtkowski, T. S. Kern, and K. Palczewski. Image registration and averaging of low laser power two-photon fluorescence images of mouse retina. *Biomedical Optics Express*, 7(7):2671, 2016.
- [4] J. Anaya and A. Barbu. Renoir-a dataset for real low-light image noise reduction. *arXiv preprint arXiv:1409.8230*, 2014.
- [5] U. Bal. Dual tree complex wavelet transform based denoising of optical microscopy images. *Biomedical Optics Express*, 3(12):3231, 2012.
- [6] Bo Zhang, J. Fadili, and J. Starck. Wavelets, Ridgelets, and Curvelets for Poisson Noise Removal. *IEEE Transactions on Image Processing*, 17(7):1093–1108, 2008.
- [7] A. Buades, B. Coll, and J.-M. Morel. A non-local algorithm for image denoising. In *Computer Vision and Pattern Recognition, 2005. CVPR 2005. IEEE Computer Society Conference on*, volume 2, pages 60–65. IEEE, 2005.
- [8] J. Chen, J. Chen, H. Chao, and M. Yang. Image blind denoising with generative adversarial network based noise modeling. In *Proceedings of the IEEE Conference on Computer Vision and Pattern Recognition*, pages 3155–3164, 2018.
- [9] K. Dabov, A. Foi, V. Katkovnik, and K. Egiazarian. Image Denoising by Sparse 3-D Transform-Domain Collaborative Filtering. *IEEE Transactions on Image Processing*, 16(8):2080–2095, 2007.
- [10] W. Denk, J. H. Strickler, and W. W. Webb. Two-photon laser scanning fluorescence microscopy. *Science*, 248(4951):73–76, 1990.
- [11] A. Foi, M. Trimeche, V. Katkovnik, and K. Egiazarian. Practical Poissonian-Gaussian noise modeling and fitting for single-image raw-data. *IEEE Transactions on Image Processing*, 17(10):1737–1754, 2008.
- [12] S. Gu, L. Zhang, W. Zuo, and X. Feng. Weighted nuclear norm minimization with application to image denoising. In *Proceedings of the IEEE Conference on Computer Vision and Pattern Recognition*, pages 2862–2869, 2014.
- [13] S. Guo, Z. Yan, K. Zhang, W. Zuo, and L. Zhang. Toward convolutional blind denoising of real photographs. *arXiv preprint arXiv:1807.04686*, 2018.
- [14] S. Lefkimmiatis. Universal denoising networks: A novel cnn-based network architecture for image denoising. *arXiv preprint arXiv:1711.07807*, 2017.
- [15] J. Lehtinen, J. Munkberg, J. Hasselgren, S. Laine, T. Karras, M. Aittala, and T. Aila. Noise2noise: Learning image restoration without clean data. *arXiv preprint arXiv:1803.04189*, 2018.
- [16] J. W. Lichtman and J.-A. Conchello. Fluorescence microscopy. *Nature Methods*, 2(12):910–919, 2005.
- [17] J. Long, E. Shelhamer, and T. Darrell. Fully convolutional networks for semantic segmentation. In *Proceedings of the IEEE Conference on Computer Vision and Pattern Recognition*, pages 3431–3440, 2015.
- [18] F. Luisier, T. Blu, and M. Unser. Image Denoising in Mixed Poisson-Gaussian Noise. *IEEE Transactions on Image Processing*, 20(3):696–708, 2011.
- [19] M. Makitalo and A. Foi. A Closed-Form Approximation of the Exact Unbiased Inverse of the Anscombe Variance-Stabilizing Transformation. *IEEE Transactions on Image Processing*, 20(9):2697–2698, 2011.
- [20] M. Makitalo and A. Foi. Optimal Inversion of the Generalized Anscombe Transformation for Poisson-Gaussian Noise. *IEEE Transactions on Image Processing*, 22(1):91–103, 2013.
- [21] W. Meinli, J.-C. Olivo-Marin, and E. D. Angelini. Denoising of Microscopy Images: A Review of the State-of-the-Art, and a New Sparsity-Based Method. *IEEE Transactions on Image Processing*, 27(8):3842–3856, 2018.
- [22] J. Pawley. *Handbook of biological confocal microscopy*. Springer Science & Business Media, 2010.
- [23] T. Plötz and S. Roth. Benchmarking denoising algorithms with real photographs. In *Computer Vision and Pattern Recognition (CVPR), 2017 IEEE Conference on*, pages 2750–2759. IEEE, 2017.
- [24] T. Plötz and S. Roth. Neural nearest neighbors networks. *arXiv preprint arXiv:1810.12575*, 2018.
- [25] P. J. Verveer, M. J. Gemkow, and T. M. Jovin. A comparison of image restoration approaches applied to three-dimensional confocal and wide-field fluorescence microscopy. *Journal of Microscopy*, 193(1):50–61, 1999.
- [26] Z. Wang, A. C. Bovik, H. R. Sheikh, and E. P. Simoncelli. Image quality assessment: from error visibility to structural similarity. *IEEE Transactions on Image Processing*, 13(4):600–612, 2004.
- [27] J. Xu, H. Li, Z. Liang, D. Zhang, and L. Zhang. Real-world noisy image denoising: A new benchmark. *arXiv preprint arXiv:1804.02603*, 2018.

- [28] K. Zhang, W. Zuo, Y. Chen, D. Meng, and L. Zhang. Beyond a gaussian denoiser: Residual learning of deep cnn for image denoising. *IEEE Transactions on Image Processing*, 26(7):3142–3155, 2017.
- [29] K. Zhang, W. Zuo, and L. Zhang. Ffdnet: Toward a fast and flexible solution for CNN based image denoising. *IEEE Transactions on Image Processing*, 2018.
- [30] Y. Zhang, G. D. Vigil, L. Cao, A. A. Khan, D. Benirschke, T. Ahmed, P. Fay, and S. S. Howard. Saturation-compensated measurements for fluorescence lifetime imaging microscopy. *Optics Letters*, 42(1):155, 2017.
- [31] D. Zoran and Y. Weiss. From learning models of natural image patches to whole image restoration. In *Computer Vision (ICCV), 2011 IEEE International Conference on*, pages 479–486. IEEE, 2011.

A. Supplementary Material

A.1. Benchmark results on separate test set

Here we show the benchmark results on the 19-th FOV (which is pre-selected as the test set) for each imaging configuration and each noise level, which contains 50 noise realizations in each case. The results are organized in Table 3 (confocal), 4 (two-photon), and 5 (wide-field). For all test cases, deep learning based denoising methods almost dominate over traditional methods.

A.2. Network architecture and training details

Network We try our best to maintain the same network structure of DnCNN and Noise2Noise as the original papers. More specifically, the only modification to Noise2Noise model is the insertion of batch normalization layer after each convolution layer and the addition of Tanh activation before the network output. For more details, please refer to the official implementations of DnCNN² and Noise2Noise³.

Training Input images are of size 256×256 , normalized to the range $[-0.5, 0.5]$. Adam optimizer is used with hyperparameters $\beta_1 = 0.9, \beta_2 = 0.99$, weight decay 0.0. The learning rate scheduling follows the one cycle policy⁴, with maximum learning rate to be 0.001, initial learning rate to be 1/10 of the maximum rate, then linearly increasing the learning rate to the maximum within 0.3 of the total epochs, then cosine annealing of the learning rate to 1/10⁵ of the maximum learning rate. The model is trained for 300 epochs. All the settings above are the same for both DnCNN and Noise2Noise.

The minibatch size is 12 for DnCNN and 16 for Noise2Noise. For DnCNN, we randomly sample 3 noisy images of size 512×512 from the training set and crop each large image into 4 non-overlapping patches of size 256×256 , thus the mini-batch size is actually 12. For Noise2Noise, we randomly select 4 pairs of large noisy images of size 512×512 where each pair (as noisy input and noisy output for training) is guaranteed belong to the same FOV. Four crop is applied as well thus the minibatch size is 16.

²<https://github.com/csxn/DnCNN>

³<https://github.com/NVlabs/noise2noise>

⁴https://github.com/fastai/fastai/blob/master/fastai/callbacks/one_cycle.py

Confocal Microscopy		Number of raw images for averaging					Time
Samples	Methods	1	2	4	8	16	
BPAE (Nuclei)	VST+NLM	37.35 / 0.9656	38.20 / 0.9730	39.31 / 0.9810	41.11 / 0.9862	43.68 / 0.9906	129.92 s
	VST+BM3D	38.45 / 0.9732	39.59 / 0.9786	40.95 / 0.9853	42.37 / 0.9889	44.37 / 0.9918	5.13 s
	VST+KSVD	38.15 / 0.9699	39.48 / 0.9773	40.92 / 0.9850	42.30 / 0.9888	44.32 / 0.9919	65.90 s
	VST+KSVD(D)	37.77 / 0.9679	39.26 / 0.9762	40.82 / 0.9846	42.28 / 0.9887	44.28 / 0.9918	20.25 s
	VST+KSVD(G)	38.07 / 0.9694	39.39 / 0.9767	40.87 / 0.9847	42.28 / 0.9886	44.20 / 0.9917	17.11 s
	VST+EPLL	38.38 / 0.9731	39.47 / 0.9785	40.85 / 0.9854	42.35 / 0.9891	44.42 / 0.9920	246.47 s
	VST+WNNM	38.43 / 0.9734	39.55 / 0.9784	40.91 / 0.9851	42.28 / 0.9885	44.21 / 0.9914	417.71 s
	PURE-LET	37.15 / 0.9583	38.55 / 0.9688	40.15 / 0.9795	41.55 / 0.9843	43.51 / 0.9887	2.43 s
	DnCNN	38.91 / 0.9795	40.23 / 0.9834	41.62 / 0.9872	43.07 / 0.9903	44.97 / 0.9930	32.48 s
	Noise2Noise	39.13 / 0.9771	40.29 / 0.9823	41.47 / 0.9858	42.73 / 0.9885	44.21 / 0.9907	9.16 s
BPAE (F-actin)	VST+NLM	32.80 / 0.8419	34.28 / 0.8893	35.76 / 0.9237	37.37 / 0.9462	39.39 / 0.9624	134.04 s
	VST+BM3D	34.07 / 0.8880	35.38 / 0.9168	36.74 / 0.9395	38.15 / 0.9556	39.80 / 0.9675	6.42 s
	VST+KSVD	33.33 / 0.8565	34.81 / 0.8985	36.25 / 0.9291	37.65 / 0.9484	39.17 / 0.9614	287.22 s
	VST+KSVD(D)	32.88 / 0.8412	34.49 / 0.8892	36.07 / 0.9245	37.55 / 0.9460	39.11 / 0.9598	64.16 s
	VST+KSVD(G)	33.08 / 0.8465	34.62 / 0.8914	36.14 / 0.9248	37.60 / 0.9457	39.17 / 0.9595	47.82 s
	VST+EPLL	34.07 / 0.8892	35.49 / 0.9207	36.94 / 0.9441	38.48 / 0.9604	40.35 / 0.9725	317.13 s
	VST+WNNM	33.94 / 0.8809	35.29 / 0.9126	36.59 / 0.9362	37.84 / 0.9515	39.21 / 0.9621	415.91 s
	PURE-LET	33.50 / 0.8776	34.75 / 0.9066	35.98 / 0.9283	37.16 / 0.9433	38.18 / 0.9505	2.66 s
	DnCNN	34.21 / 0.9029	35.62 / 0.9311	37.07 / 0.9512	38.66 / 0.9665	40.75 / 0.9791	33.81 s
	Noise2Noise	34.33 / 0.9025	35.63 / 0.9289	36.92 / 0.9480	38.30 / 0.9625	39.92 / 0.9736	9.30 s
BPAE (Mito)	VST+NLM	35.79 / 0.9279	37.27 / 0.9518	38.93 / 0.9673	40.89 / 0.9781	43.36 / 0.9865	130.14 s
	VST+BM3D	37.43 / 0.9489	38.82 / 0.9632	40.27 / 0.9742	41.80 / 0.9817	43.78 / 0.9879	5.92 s
	VST+KSVD	36.97 / 0.9378	38.49 / 0.9575	39.98 / 0.9712	41.48 / 0.9798	43.33 / 0.9865	241.33 s
	VST+KSVD(D)	36.55 / 0.9305	38.25 / 0.9537	39.89 / 0.9695	41.50 / 0.9792	43.42 / 0.9864	60.91 s
	VST+KSVD(G)	36.93 / 0.9368	38.59 / 0.9579	40.18 / 0.9720	41.71 / 0.9806	43.59 / 0.9871	42.51 s
	VST+EPLL	37.56 / 0.9515	38.95 / 0.9653	40.41 / 0.9757	41.94 / 0.9828	43.98 / 0.9887	312.86 s
	VST+WNNM	37.46 / 0.9486	38.91 / 0.9638	40.34 / 0.9745	41.80 / 0.9816	43.67 / 0.9875	502.87 s
	PURE-LET	36.87 / 0.9433	38.12 / 0.9568	39.47 / 0.9678	40.95 / 0.9764	42.73 / 0.9834	2.70 s
	DnCNN	37.89 / 0.9586	39.30 / 0.9702	40.68 / 0.9781	42.14 / 0.9841	44.00 / 0.9894	32.08 s
	Noise2Noise	37.74 / 0.9549	39.13 / 0.9675	40.47 / 0.9756	41.78 / 0.9813	43.22 / 0.9859	8.77 s
Zebrafish Embryo	VST+NLM	28.23 / 0.7895	31.47 / 0.8593	34.00 / 0.9078	35.72 / 0.9328	37.58 / 0.9482	145.64 s
	VST+BM3D	32.00 / 0.8854	33.75 / 0.9102	35.30 / 0.9301	36.78 / 0.9443	38.32 / 0.9546	6.29 s
	VST+KSVD	29.04 / 0.8203	32.17 / 0.8740	34.58 / 0.9167	36.31 / 0.9388	37.86 / 0.9519	60.01 s
	VST+KSVD(D)	28.87 / 0.8184	31.42 / 0.8647	33.97 / 0.9093	35.97 / 0.9350	37.74 / 0.9504	12.54 s
	VST+KSVD(G)	29.03 / 0.8201	31.88 / 0.8701	34.34 / 0.9133	36.26 / 0.9374	38.04 / 0.9520	9.93 s
	VST+EPLL	31.62 / 0.8678	33.66 / 0.9048	35.34 / 0.9298	36.92 / 0.9460	38.61 / 0.9574	317.67 s
	VST+WNNM	30.94 / 0.8654	33.43 / 0.9048	35.23 / 0.9284	36.74 / 0.9432	38.14 / 0.9527	615.40 s
	PURE-LET	30.03 / 0.8019	32.48 / 0.8817	33.84 / 0.8960	35.65 / 0.9254	37.15 / 0.9394	2.59 s
	DnCNN	32.44 / 0.9025	34.16 / 0.9267	35.75 / 0.9425	37.28 / 0.9548	39.07 / 0.9659	32.25 s
	Noise2Noise	32.93 / 0.9076	34.37 / 0.9270	35.71 / 0.9410	37.06 / 0.9523	38.65 / 0.9625	9.44 s
Mouse Brain	VST+NLM	36.31 / 0.9534	37.53 / 0.9632	38.95 / 0.9706	40.87 / 0.9763	43.37 / 0.9819	131.08 s
	VST+BM3D	37.95 / 0.9637	39.47 / 0.9704	41.09 / 0.9765	42.73 / 0.9811	44.52 / 0.9847	6.24 s
	VST+KSVD	37.46 / 0.9587	39.24 / 0.9684	40.94 / 0.9757	42.55 / 0.9807	44.24 / 0.9846	85.33 s
	VST+KSVD(D)	36.67 / 0.9544	38.68 / 0.9659	40.63 / 0.9746	42.43 / 0.9804	44.26 / 0.9846	21.95 s
	VST+KSVD(G)	37.30 / 0.9582	39.15 / 0.9681	40.93 / 0.9757	42.65 / 0.9808	44.49 / 0.9849	17.89 s
	VST+EPLL	37.92 / 0.9640	39.50 / 0.9710	41.18 / 0.9772	42.87 / 0.9818	44.73 / 0.9855	320.98 s
	VST+WNNM	37.86 / 0.9624	39.47 / 0.9698	41.08 / 0.9761	42.62 / 0.9804	44.17 / 0.9837	456.09 s
	PURE-LET	36.60 / 0.9359	38.10 / 0.9477	40.06 / 0.9650	41.75 / 0.9739	43.29 / 0.9791	2.54 s
	DnCNN	38.15 / 0.9672	39.78 / 0.9741	41.41 / 0.9794	43.11 / 0.9841	45.20 / 0.9887	33.86 s
	Noise2Noise	38.19 / 0.9665	39.77 / 0.9735	41.28 / 0.9787	42.83 / 0.9831	44.56 / 0.9869	8.69 s

Table 3. Denoising performance of confocal microscopy images (the 19-th FOV of each imaging configuration). PSNR (dB), SSIM, and denoising time (seconds) are obtained by averaging over 50 noise realizations through imaging experiments.

Two-Photon Microscopy		Number of raw images for averaging					Time
Samples	Methods	1	2	4	8	16	
BPAE (Nuclei)	VST+NLM	31.34 / 0.9173	32.13 / 0.9286	32.95 / 0.9390	34.14 / 0.9482	37.35 / 0.9571	137.27 s
	VST+BM3D	32.02 / 0.9297	32.70 / 0.9382	33.43 / 0.9458	34.60 / 0.9526	37.77 / 0.9592	5.58 s
	VST+KSVD	31.71 / 0.9227	32.55 / 0.9352	33.37 / 0.9453	34.55 / 0.9535	37.70 / 0.9613	42.51 s
	VST+KSVD(D)	31.48 / 0.9195	32.33 / 0.9323	33.23 / 0.9438	34.48 / 0.9529	37.69 / 0.9612	10.77 s
	VST+KSVD(G)	31.70 / 0.9225	32.52 / 0.9347	33.34 / 0.9448	34.55 / 0.9533	37.75 / 0.9613	8.12 s
	VST+EPLL	32.00 / 0.9313	32.70 / 0.9404	33.48 / 0.9483	34.69 / 0.9552	37.95 / 0.9618	284.32 s
	VST+WNNM	32.01 / 0.9298	32.68 / 0.9383	33.41 / 0.9460	34.55 / 0.9524	37.62 / 0.9585	487.02 s
	PURE-LET	31.62 / 0.9101	32.27 / 0.9198	32.88 / 0.9231	33.97 / 0.9312	36.92 / 0.9439	2.68 s
	DnCNN	31.59 / 0.9250	32.46 / 0.9421	33.38 / 0.9513	34.75 / 0.9598	38.30 / 0.9705	32.42 s
	Noise2Noise	32.44 / 0.9354	33.21 / 0.9434	34.04 / 0.9509	35.19 / 0.9590	38.22 / 0.9685	8.14 s
BPAE (F-actin)	VST+NLM	30.26 / 0.7176	31.43 / 0.7799	32.70 / 0.8404	34.24 / 0.8912	37.04 / 0.9297	229.93 s
	VST+BM3D	31.59 / 0.8037	32.52 / 0.8442	33.56 / 0.8813	34.91 / 0.9139	37.56 / 0.9408	5.89 s
	VST+KSVD	30.67 / 0.7381	31.84 / 0.7992	33.10 / 0.8560	34.54 / 0.8995	37.07 / 0.9304	163.48 s
	VST+KSVD(D)	30.43 / 0.7261	31.52 / 0.7833	32.83 / 0.8438	34.38 / 0.8936	37.00 / 0.9279	30.14 s
	VST+KSVD(G)	30.57 / 0.7325	31.69 / 0.7904	32.97 / 0.8485	34.48 / 0.8952	37.09 / 0.9284	24.08 s
	VST+EPLL	31.48 / 0.7950	32.56 / 0.8456	33.72 / 0.8889	35.19 / 0.9237	38.09 / 0.9507	287.27 s
	VST+WNNM	31.24 / 0.7778	32.30 / 0.8278	33.41 / 0.8723	34.76 / 0.9082	37.25 / 0.9345	506.98 s
	PURE-LET	31.19 / 0.7858	32.09 / 0.8267	33.19 / 0.8705	34.53 / 0.9055	36.85 / 0.9295	2.62 s
	DnCNN	31.52 / 0.8222	32.67 / 0.8685	33.92 / 0.9059	35.47 / 0.9368	38.68 / 0.9643	31.65 s
	Noise2Noise	32.00 / 0.8257	33.10 / 0.8701	34.19 / 0.9048	35.59 / 0.9342	38.46 / 0.9596	8.53 s
BPAE (Mito)	VST+NLM	35.11 / 0.8525	36.73 / 0.8917	38.66 / 0.9290	40.68 / 0.9554	43.49 / 0.9738	208.28 s
	VST+BM3D	37.52 / 0.9130	38.72 / 0.9338	40.09 / 0.9511	41.62 / 0.9648	43.97 / 0.9766	5.49 s
	VST+KSVD	35.75 / 0.8679	37.34 / 0.9039	39.21 / 0.9367	40.98 / 0.9576	43.29 / 0.9725	97.25 s
	VST+KSVD(D)	35.61 / 0.8648	36.96 / 0.8961	38.77 / 0.9295	40.66 / 0.9536	43.12 / 0.9710	19.32 s
	VST+KSVD(G)	35.74 / 0.8675	37.25 / 0.9019	39.16 / 0.9354	41.07 / 0.9579	43.57 / 0.9737	14.39 s
	VST+EPLL	37.29 / 0.9065	38.81 / 0.9348	40.38 / 0.9549	42.05 / 0.9689	44.58 / 0.9800	291.54 s
	VST+WNNM	36.68 / 0.8929	38.30 / 0.9250	39.90 / 0.9481	41.51 / 0.9636	43.77 / 0.9754	525.45 s
	PURE-LET	36.88 / 0.8946	38.01 / 0.9179	38.70 / 0.9276	40.12 / 0.9459	42.27 / 0.9637	2.77 s
	DnCNN	38.15 / 0.9251	39.46 / 0.9460	40.87 / 0.9616	42.51 / 0.9738	45.32 / 0.9845	33.85 s
	Noise2Noise	38.11 / 0.9241	39.38 / 0.9450	40.77 / 0.9606	42.37 / 0.9727	44.82 / 0.9825	8.72 s
Mouse Brain	VST+NLM	32.80 / 0.9134	33.88 / 0.9237	34.88 / 0.9317	36.31 / 0.9384	38.96 / 0.9449	211.65 s
	VST+BM3D	33.81 / 0.9246	34.78 / 0.9317	35.77 / 0.9379	36.97 / 0.9431	39.39 / 0.9481	6.14 s
	VST+KSVD	33.35 / 0.9183	34.47 / 0.9288	35.60 / 0.9374	36.85 / 0.9442	39.27 / 0.9509	79.00 s
	VST+KSVD(D)	32.89 / 0.9147	34.14 / 0.9264	35.43 / 0.9362	36.79 / 0.9437	39.26 / 0.9507	13.64 s
	VST+KSVD(G)	33.34 / 0.9179	34.50 / 0.9285	35.66 / 0.9372	36.94 / 0.9441	39.42 / 0.9508	9.83 s
	VST+EPLL	33.86 / 0.9262	34.86 / 0.9339	35.86 / 0.9403	37.11 / 0.9456	39.61 / 0.9506	286.50 s
	VST+WNNM	33.79 / 0.9254	34.75 / 0.9323	35.74 / 0.9386	36.91 / 0.9435	39.22 / 0.9480	512.61 s
	PURE-LET	32.86 / 0.8812	33.47 / 0.8720	34.42 / 0.8769	35.49 / 0.8878	37.40 / 0.8997	2.84 s
	DnCNN	33.67 / 0.9068	34.95 / 0.9290	36.10 / 0.9413	37.43 / 0.9507	40.30 / 0.9630	34.49 s
	Noise2Noise	34.33 / 0.9249	35.32 / 0.9335	36.25 / 0.9410	37.46 / 0.9499	39.89 / 0.9609	8.49 s

Table 4. Denoising performance of two-photon microscopy images (the 19-th FOV of each imaging configuration). PSNR (dB), SSIM, and denoising time (seconds) are obtained by averaging over 50 noise realizations through imaging experiments.

Wide-Field Microscopy		Number of raw images for averaging					Time
Samples	Methods	1	2	4	8	16	
BPAE (Nuclei)	VST+NLM	25.53 / 0.3875	28.49 / 0.5548	31.36 / 0.7122	34.33 / 0.8397	37.74 / 0.9264	138.54 s
	VST+BM3D	26.22 / 0.4339	29.16 / 0.6020	31.99 / 0.7511	34.91 / 0.8650	38.25 / 0.9386	6.13 s
	VST+KSVD	26.38 / 0.4459	29.31 / 0.6132	32.10 / 0.7577	34.99 / 0.8681	38.30 / 0.9397	1348.61 s
	VST+KSVD(D)	26.41 / 0.4489	29.33 / 0.6152	32.11 / 0.7590	35.00 / 0.8688	38.30 / 0.9398	183.82 s
	VST+KSVD(G)	26.40 / 0.4533	29.32 / 0.6182	32.10 / 0.7604	34.98 / 0.8693	38.28 / 0.9399	170.94 s
	VST+EPLL	26.06 / 0.4244	29.00 / 0.5923	31.86 / 0.7440	34.79 / 0.8601	38.15 / 0.9365	354.13 s
	VST+WNNM	26.36 / 0.4440	29.29 / 0.6116	32.11 / 0.7581	35.01 / 0.8690	38.32 / 0.9402	420.74 s
	PURE-LET	26.13 / 0.4258	29.05 / 0.5931	31.89 / 0.7442	34.79 / 0.8593	38.07 / 0.9341	2.49 s
	DnCNN	32.94 / 0.8665	35.42 / 0.9225	37.04 / 0.9440	38.30 / 0.9539	40.07 / 0.9645	30.82 s
	Noise2Noise	34.19 / 0.9047	36.13 / 0.9339	37.42 / 0.9475	38.58 / 0.9563	40.12 / 0.9650	9.03 s
BPAE (F-actin)	VST+NLM	23.93 / 0.3370	27.02 / 0.4988	30.21 / 0.6672	33.58 / 0.8096	37.67 / 0.9150	132.00 s
	VST+BM3D	24.72 / 0.3792	27.84 / 0.5467	31.02 / 0.7084	34.36 / 0.8367	38.27 / 0.9258	5.66 s
	VST+KSVD	24.94 / 0.3910	28.03 / 0.5575	31.22 / 0.7178	34.54 / 0.8426	38.48 / -0.9292	1343.88 s
	VST+KSVD(D)	25.01 / 0.3965	28.11 / 0.5629	31.28 / 0.7213	34.59 / 0.8445	38.51 / 0.9297	175.55 s
	VST+KSVD(G)	25.04 / 0.4036	28.13 / 0.5683	31.30 / 0.7245	34.60 / 0.8458	38.50 / 0.9299	156.79 s
	VST+EPLL	24.55 / 0.3711	27.70 / 0.5393	30.88 / 0.7018	34.24 / 0.8331	38.16 / 0.9241	352.19 s
	VST+WNNM	24.94 / 0.3900	28.01 / 0.5560	31.17 / 0.7154	34.48 / 0.8406	38.36 / 0.9272	438.09 s
	PURE-LET	24.67 / 0.3736	27.75 / 0.5393	30.90 / 0.7012	34.18 / 0.8306	37.64 / 0.9134	2.49 s
	DnCNN	32.49 / 0.8044	34.49 / 0.8561	36.00 / 0.8864	37.60 / 0.9151	39.47 / 0.9408	31.25 s
	Noise2Noise	33.04 / 0.8200	34.53 / 0.8548	36.03 / 0.8872	37.70 / 0.9182	39.64 / 0.9451	9.15 s
BPAE (Mito)	VST+NLM	26.26 / 0.4134	29.35 / 0.5850	32.55 / 0.7418	35.96 / 0.8610	39.93 / 0.9389	134.42 s
	VST+BM3D	26.93 / 0.4611	30.03 / 0.6312	33.24 / 0.7778	36.65 / 0.8831	40.58 / 0.9487	5.97 s
	VST+KSVD	27.11 / 0.4737	30.20 / 0.6417	33.38 / 0.7845	36.76 / 0.8863	40.70 / 0.9504	1247.01 s
	VST+KSVD(D)	27.14 / 0.4768	30.22 / 0.6440	33.40 / 0.7859	36.78 / 0.8869	40.69 / 0.9504	172.92 s
	VST+KSVD(G)	27.13 / 0.4804	30.22 / 0.6464	33.40 / 0.7870	36.76 / 0.8872	40.66 / 0.9503	161.30 s
	VST+EPLL	26.80 / 0.4524	29.91 / 0.6233	33.12 / 0.7721	36.51 / 0.8791	40.46 / 0.9471	345.91 s
	VST+WNNM	27.08 / 0.4709	30.17 / 0.6400	33.37 / 0.7841	36.77 / 0.8866	40.69 / 0.9502	430.37 s
	PURE-LET	26.85 / 0.4528	29.94 / 0.6231	33.13 / 0.7709	36.49 / 0.8777	40.27 / 0.9440	2.56 s
	DnCNN	34.54 / 0.8881	37.05 / 0.9253	38.93 / 0.9429	40.70 / 0.9557	42.72 / 0.9679	32.58 s
	Noise2Noise	34.71 / 0.8926	36.83 / 0.9215	38.80 / 0.9408	40.66 / 0.9552	42.62 / 0.9677	9.51 s

Table 5. Denoising performance of wide-field microscopy images (the 19-th FOV of each imaging configuration). PSNR (dB), SSIM, and denoising time (seconds) are obtained by averaging over 50 noise realizations through imaging experiments.