

Title: Methods to simplify cooling of liquid Helium cryostats

Authors: Rafael Álvarez Montoya¹, Sara Delgado¹, José Castilla¹, José Navarrete², Nuria Díaz Contreras², Juan Ramón Marijuan², Víctor Barrena,¹ Isabel Guillamón¹ and Hermann Suderow¹

Affiliations: ¹ Laboratorio de Bajas Temperaturas, Departamento de Física de la Materia Condensada

Instituto de Ciencia de Materiales Nicolás Cabrera and Condensed Matter Physics Center (IFIMAC), Facultad de Ciencias

Universidad Autónoma de Madrid, 28049 Madrid, Spain

² Segainvex

Universidad Autónoma de Madrid, 28049 Madrid, Spain

Contact email: hermann.suderow@uam.es

Abstract: Liquid Helium is used widely, from hospitals to characterization of materials at low temperatures. Many experiments at low temperatures require liquid Helium, particularly when vibration isolation precludes the use of cryocoolers and when one needs to cool heavy equipment such as superconducting coils. Here we describe methods to simplify the operations required to use liquid Helium by eliminating the use of high pressure bottles, avoiding blockage and improving heating and cooling rates. First we show a simple and very low cost method to transfer liquid Helium from a transport container into a cryostat that uses a manual pump having pumping and pressurizing ports, giving a liquid Helium transfer rate of about 100 liters an hour. Second, we describe a closed cycle circuit of Helium gas cooled in an external liquid nitrogen bath that allows precooling a cryogenic experiment without inserting liquid nitrogen into the cryostat, eliminating problems associated to the presence of nitrogen around superconducting magnets. And third, we show a sliding seal assembly and an inner vacuum chamber design that allows inserting large experiments into liquid Helium.

Keywords: Liquid Helium cryogenics. Superconducting solenoids. Cryogenic scanning probe microscopy. Materials characterization at low temperatures and high magnetic fields.

Specifications tables:

Hardware name	Helium transfer system based on manual pump
Subject area	<ul style="list-style-type: none"> • Materials Science • Physical property measurements • Cryogenics of magnetic resonance equipment
Hardware type	<ul style="list-style-type: none"> • Device to transfer liquid Helium • Maintenance and preparation of systems for physical property measurements at low temperatures
Open source license	CERN OHL
Cost of hardware	30 \$
Source file repository	https://osf.io/e6k7r/

Hardware name	Closed cycle Helium gas precooling system
Subject area	<ul style="list-style-type: none"> • Materials Science • Physical property measurements
Hardware type	<ul style="list-style-type: none"> • Precooling equipment to liquid nitrogen temperatures • Physical property measurements at low temperatures
Open source license	CERN OHL
Cost of hardware	1500 \$
Source file repository	https://osf.io/e6k7r/

Hardware name	Sliding seal and inner vacuum chamber for rapid sample exchange
Subject area	<ul style="list-style-type: none"> • Materials Science • Physical property measurements
Hardware type	<ul style="list-style-type: none"> • Precooling equipment to liquid nitrogen temperatures • Physical property measurements at low temperatures
Open source license	<i>CERN OHL</i>
Cost of hardware	200 \$
Source file repository	https://osf.io/e6k7r/

1. Hardware in context

Much of the innovative work in cryogenics is now devoted into dry or closed-cycle cryogenic systems, driven by the development of thermoacoustic systems and the improvement of thermodynamic machines [1]. However, the improvements have not displaced liquid Helium, but considerably facilitated its development and use in different fields. For example, magnetic resonance imaging cryostats include coolers providing improved shielding that nearly zero the consumption of liquid Helium[2, 3]. In research laboratories, large systems including superconducting magnets require most often liquid Helium, because cooling using the huge enthalphy of Helium gas is unmatched with respect to any other cooling method. In addition, there is an increasing need for experiments in which the mechanical vibrations are held to a minimum, for which the use of liquid Helium is unavoidable[4, 5, 6, 7, 8, 9, 10, 11]. Therefore, research in wet cryogenics is needed to improve the performance of systems using liquid Helium.

In particular, cooling of large cryostats bears some inherent problems related to the manipulation of liquid Helium. There are many laboratory books and notes that are there to help the users, in addition to usual manuals provided by companies delivering cryostats [12, 13]. Most of these manuals explain transfer of liquid Helium and sample turn-around methods. But several problems are repeatedly mentioned. These are about handling of liquid Helium bottles and transfer methods, inserting and cooling to liquid Helium temperatures large inserts, often containing a dilution refrigerator unit, and precooling to liquid nitrogen temperatures. For example, as we explain below, tubes can be blocked when liquid nitrogen has been used and not fully removed before transferring liquid Helium, or Helium leaks appear when inserting large systems into liquid Helium.

Here we discuss possible solutions to these issues, which we developed in the use of a set of five large dilution refrigerator units that are equipped with Scanning Tunneling Microscopes. The microscopes have been partially described in detail in Ref.[14]. A cryogenic set-up with a three axis vector magnet and associated electronics for the magnet and the microscope has been described in Ref.[15]. The solutions described here have been developed more recently. The Helium transfer manual pump has not been reported till now, to our knowledge. The liquid nitrogen precooling device has been reported in Ref.[16] in a realization that includes a specific compressor and heat exchanging system. Here we show a much simpler realization, using commercially available

equipment. The improvements in sliding seal assembly and vacuum chamber have not been reported to our knowledge.

2. Hardware description

Potential advantages of using the described hardware:

- Eliminating high pressure Helium equipment from cryogenic laboratories.
- Avoiding blockage and problems with transfer of liquid nitrogen into cryostats.
- Improving sample turn around time in cryogenic devices.

2.1 Helium transfer system based on manual pump

Liquid Helium is usually transferred from a transport Dewar to a cryostat using a transfer tube and applying an overpressure to the transport Dewar. An acceptable flow rate is usually of about 50 liters of liquid Helium for half an hour, which makes a volumetric flow rate $Q \approx 30$ ml/s. With a cross section of a pipe which typically is of about $d = 10\text{mm}^2$, this gives a velocity of liquid Helium close to the speed of sound. The force needed to produce Helium flow is the force needed to overcome height differences in the liquid Helium levels of both containers and to overcome friction. An overpressure of 0.2-0.4 bar in a container of typically a quarter of a square meter cross section is enough to transfer liquid Helium into the cryostat. Although the Reynolds number can be as high as 10^6 , losses due to friction do not hamper significantly the transfer process[17]. The high Reynolds number makes liquid Helium flow a rather unusual process, which is however routinely carried out in laboratories.

To produce the overpressure of, say, 0.5 bar, one can simply introduce Helium gas into the transport Dewar by some means. A few simple calculations show that the evaporation produced in the liquid by the arrival of a small amount of hot gas is by far enough to produce an overpressure in the transport Dewar. For example, in a usual 100 liters transport Dewar, there is somewhat less than a m^3 room for the gas. The latent heat of vaporization of liquid Helium is of about 20kJ/kg. The enthalpy of Helium gas between room temperature and liquid Helium temperatures is of about 1.5 MJ/kg. Cooling of a kg of gas to liquid Helium temperatures requires evaporating about 75 kg of liquid. Taking into account the different densities, this implies that to achieve the required overpressure, one just needs to evaporate about a liter of liquid and that this can be done by cooling to liquid Helium of just about 12 liters of gas. Therefore, one needs to insert into the transport Dewar a small amount of Helium gas at room temperature to produce the required overpressure.

A usual solution is a small laboratory membrane compressor, which allows in principle for automated operation, but it is often not very practical. The usual flow rate of a small devices is of about 10 liters per minute or less. In practice, the compressor is often stuck by impurities or debris from the recovery line and one needs to operate the compressor for several minutes to reach the required pressure. This implies that the gas has enough time to thermalize inside the transport Dewar and remain at the upper level, without really evaporating the liquid. Furthermore, care has to be taken to make sure that Helium is not pumped out of the recovery line when it is not needed, which would imply losses of Helium gas. Therefore, most laboratories use insted Helium gas bottles pressurized to 200 bar with a pressure reducing manometer. However, high pressure gas bottles requires significant security management and it is quite easy to loose large amounts of gas by bad manipulation.

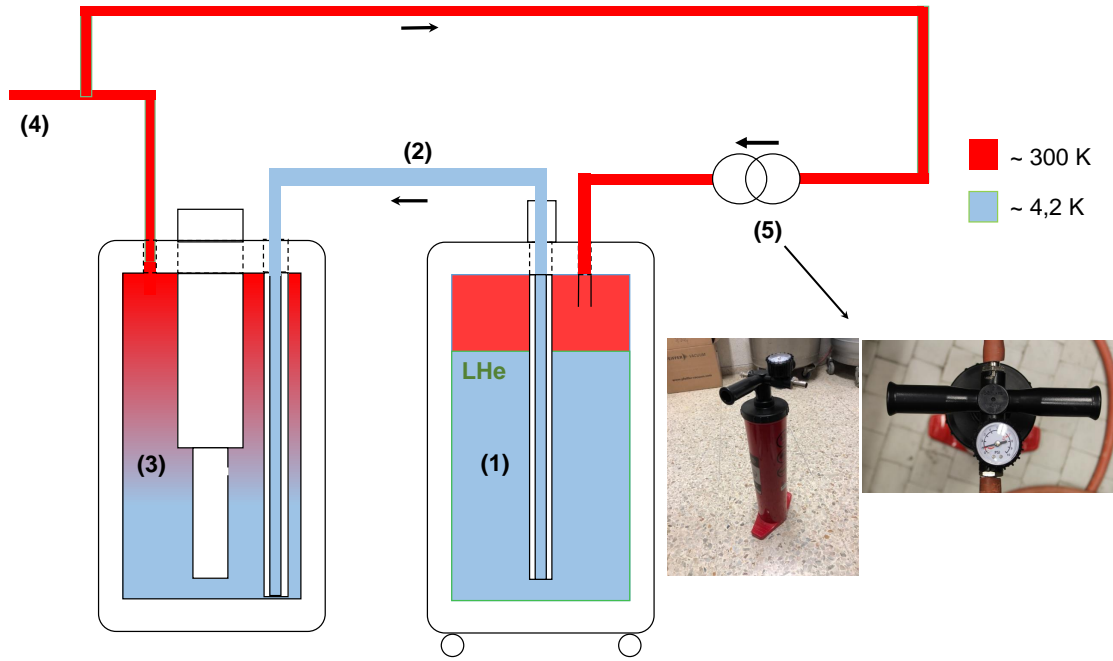


Figure 1: Scheme of the liquid Helium transfer system. We show schematically a liquid Helium transport Dewar (1), a transfer tube (2) and a cryostat (3). The cryostat is connected to the Helium recovery line (4). The transport Dewar (1) is connected to the recovery line (4) through a manual pump (5). The manual pump (5) is shown as a photo in the lower right part of the figure. The manual gas pump has clear entry and exit ports. Connectors for rubber tubes are glued to the ports. About five strokes of the manual pump are needed to fastly send Helium gas into the Dewar. The Helium gas is cooled by the liquid Helium, which evaporates and increases the pressure in the transport Dewar. The evaporation produces a small overpressure of about 0.5 bar, with which liquid Helium is transferred to the cryostat. Using this pump, 100 liters can be transferred in about an hour.

Rubber bladders can be used when slow filling rates are required[13]. The bladder is inflated with the small overpressure of the Helium Dewar. By pushing on the bladder, a few liters of Helium gas are moved inside the Dewar and the obtained flow of gas evaporates a small amount of liquid, which leads to an increase in pressure. However, the bladder is inflated and it becomes very difficult to manipulate. Small pressures can be nicely regulated this way, but when high transfer rates are needed, the bladder becomes too big to operate.

Recently, manual gas pumps with an entry and an exit port have been made available in the market. These are used to fill and empty large containers with air, such as boats or camping appliances. They pump gas volumes between 5 and 10 liters per stroke and a maximal pressure difference of order of 0.5 bar (see for example https://www.decathlon.co.uk/52-1-hand-pump-id_8243066.html). We have adapted one of such pumps to collect Helium from the recovery of the cryostat and insert the gas into the transfer Dewar. The pump has the speed and size to move enough Helium gas inside the transport Dewar and produce the required liquid Helium evaporation. We show schematically in Fig.1 the whole arrangement. A few strokes to the pump are enough to produce the required increase in the pressure of the transport Dewar and fill a cryostat to with a velocity of about 100 liters an hour.

2.2 Closed cycle Helium gas precooling system

As mentioned above, the latent heat of vaporization of liquid Helium is small (21 kJ/kg), so that precooling the system to liquid nitrogen allows reducing considerably the Helium boil-off. On the other hand, the melting point of liquid nitrogen is just at 63 K, quite close to 77 K, and the latent heat of fusion is of about 26 kJ/kg. This implies that, if the liquid nitrogen is not fully removed when the system has cooled, a sizeable amount of liquid Helium is needed to solidify the liquid nitrogen and cool down to liquid Helium temperatures. For example, if just an amount of liquid nitrogen of order of a cup, or of about 100 g of liquid nitrogen, is left in the system, 2.6 kJ are needed to solidify it, which implies evaporating about a kg, or about 10 l, of liquid Helium.

Removing all liquid nitrogen is not easy. One can use a tube inserted down to the bottom of the cryostat and an overpressure, so that nitrogen is transferred back outside the cryostat. However, the tube must be designed in such a way as to reach the bottom of the cryostat so that all nitrogen can be removed. In addition, the cryostat must be flushed and pumped several times[13]. There is always a risk of leaving small amounts of nitrogen in thin tubes, and this can lead to blockage when cooling to liquid Helium. Furthermore, nitrogen might get trapped inside superconducting magnets, which produces strain that might lead to quench during operation.

To avoid these difficulties, we have designed a precooling system as schematically described in Fig.2. Helium gas is circulated within a close circuit through a heat exchanger in liquid nitrogen, then inserted into the cryostat and flows again into the circuit after cooling the cryostat. We force Helium gas at room temperature through a serpentine immersed in liquid nitrogen and insert the outgoing cold gas using a usual Helium transfer tube into a cryostat with a superconducting coil inside. The warm Helium gas exiting the cryostat is again inserted into the circuit.

The same method has been proposed in the past[16], with, however, a rather complex heat exchanger and a multiple stage blade compressor.

To discuss the working principle of the system, let us consider an example with a system whose cooling would be equivalent to cool about 50 kg of copper. The heat of vaporization of liquid nitrogen is of about 200 kJ/kg. Density of liquid nitrogen is of about 800 kg/m³, so the vaporization of a m³ of liquid provides 160MJ. Cooling 50 kg of copper from room temperature down to liquid nitrogen requires about 3 MJ (63 kJ per kg of copper), or 20 liters of liquid nitrogen.

On the other hand, Helium gas has a specific heat at constant volume of about 3kJ/kgK. Taking a density of 0.17 kg/m³, this gives 0.5 kJ/m³K. Assuming that Helium follows the ideal gas law and the specific heat is constant with temperature, we estimate that cooling Helium gas from room temperature to liquid nitrogen requires removing 100 kJ per m³ of gas. If we use liquid nitrogen to cool the gas, we need about 0.6 liters of liquid nitrogen to cool a m³ of Helium gas to liquid nitrogen temperatures.

To cool the amount of copper mentioned above (50 kg) to liquid nitrogen, we need to evaporate at least 20 liters of liquid nitrogen. Here we do this by using Helium gas that has been previously thermalized in the liquid nitrogen. This implies that we need to circulate about 30 m³ of Helium gas through liquid nitrogen. Of course, best is to use a closed-cycle circuit. With a circulation of 10m³/h of Helium gas through the liquid nitrogen, the required cooling is achieved within some hours. This circulation rate can be readily achieved using a medium size commercial membrane compressor.

It is relatively easy to build a heat exchanger (Fig.3(a)) that allows cooling such an amount of Helium gas with liquid nitrogen. A one-dimensional heat exchanger consisting of a 1 m long tube copper tube with a diameter of the order of a cm suffices to achieve the desired cooling.

The circuit requires an expansion Helium gas by the pressure difference which is established between the input in the cryostat and the exit. The Joule-Kelvin coefficient of Helium gas is small,

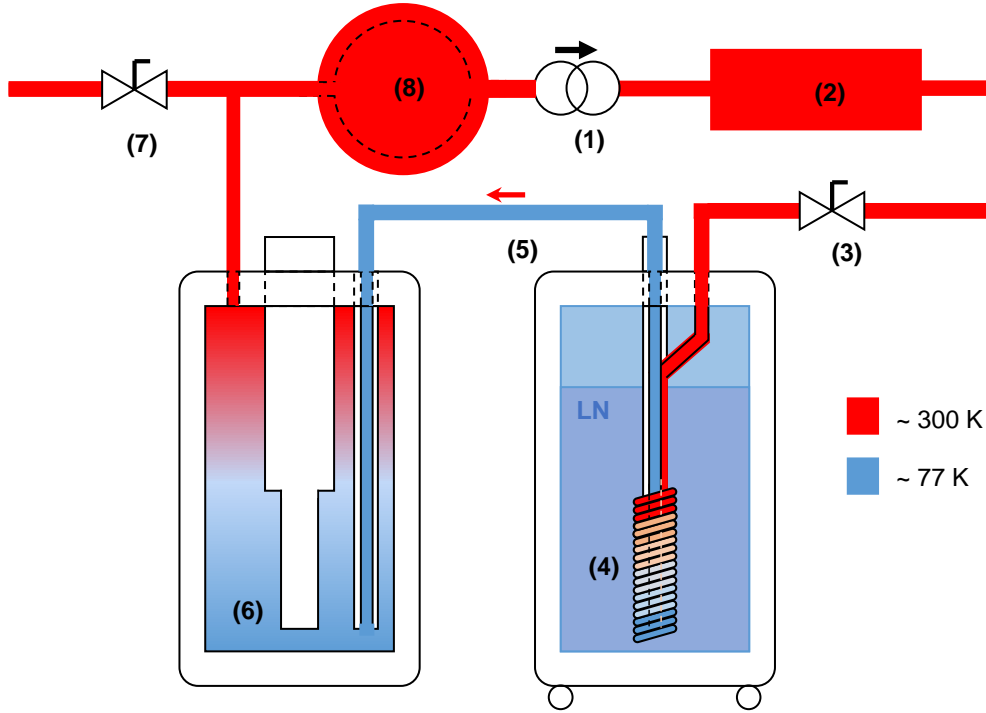


Figure 2: Scheme of the precooling system. A compressor (1) is connected to a pressure container (2). A tube connects the container (2), through a regulating valve (3) to a heat exchanger immersed into liquid nitrogen (4). Helium flows out of the heat exchanger through a usual Helium transfer tube (5) into the cryostat (6). A valve (7) connects the recovery of the cryostat to the Helium recovery line of the laboratory. This valve is closed during operation and can be used to fill the whole circuit with Helium gas. The recovery of the cryostat is connected to a balloon (8) that reduces the time needed to operate the compressor (1) and allows for continuous flow.

of about 0.06 K/bar[18], so that the expansion implies a minimal heating (about 60 mK).

We have used this system to cool a cryostat with a superconducting coil and a dilution refrigerator and could reach temperatures of about 100 K overnight (see validation results in Fig.5) and then immediately started transferring Helium.

An relevant aspect is that, if the compressor is operated all the time, it might burn down the motor of a simple commercial device. Therefore, the circuit is not a simple flow through a compressor on a closed-cycle (Fig.3(b)). The pressure difference established by the flow of Helium gas occurs at the transfer tube (Fig.2) mostly. Using Hagen-Poiseuille, we obtain about 1 bar considering a transfer tube of 2 mm diameter and 3 m length. To avoid using the compressor all times, we insert a pressurized container with a valve limiting the flow rate (Fig.3 and (2) at Fig.2). The compressor is activated whenever the pressure reaches 1.5 bar and stopped when it reaches 4 bar at the pressurized container. A valve at the exit of the pressurized container limits the pressure on the transfer line so that it remains below about 2 bar ((3) at Fig.2). In order to have enough gas to establish the overpressure each time the compressor is switched on again, we require a Helium

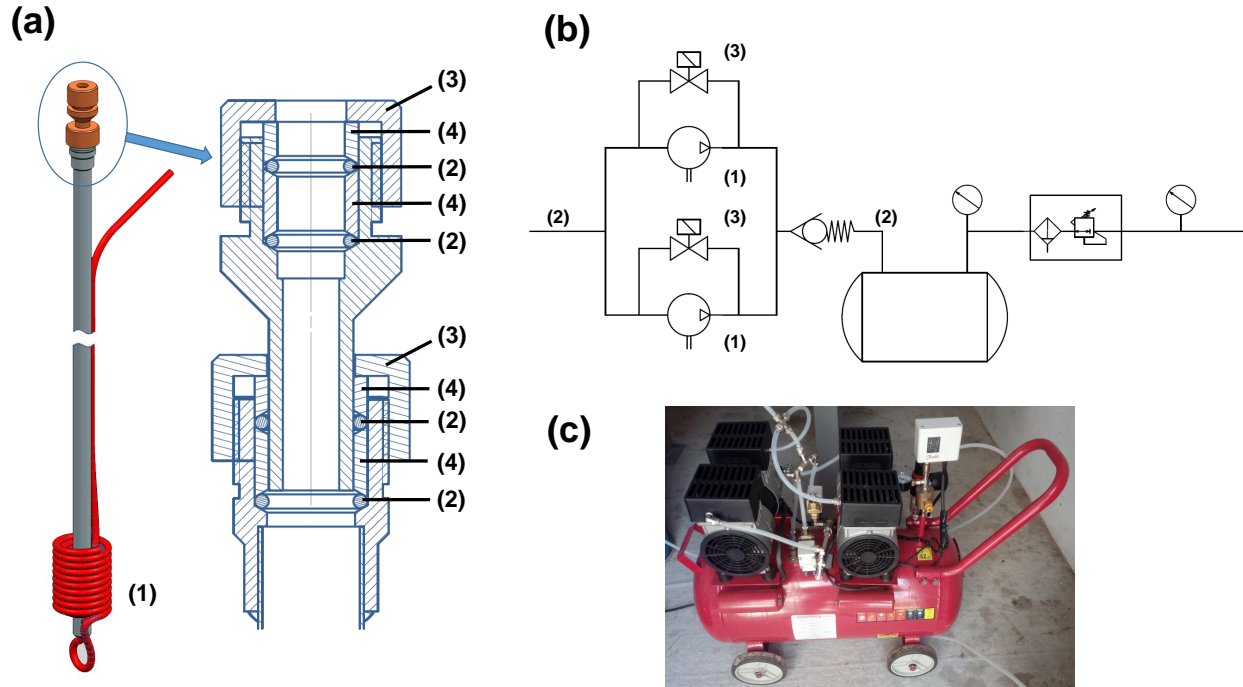


Figure 3: (a) Scheme of the heat exchanger. The heat exchanger consists of a long copper tube wound at the bottom (1) in such a way that the external diameter enters a small liquid nitrogen container. The top of the heat exchanger (the exit) has a O-ring system that allows inserting transfer tubes of the two most common diameters (9 mm and 12 mm). It consists of four O-rings (2) that are held by threaded brass caps (3) and internal brass caps (4) that are machined at an angle at the ends to compress the O-ring to the transfer tube that is to be inserted. (b) Scheme of the compressor. The system (1) is a commercially available set of two membrane compressors mounted in parallel. The entry and exit ports (2) have been sealed using Araldite. The valves (3) have been exchanged by leak-tight valves and we have added a system that allows equilibrating pressures before re-starting the compressor. (c) Photography of the modified compressor with pressure container.

balloon on the low pressure side ((8) at Fig.2). The balloon should be at least four times the volume of the pressurized container.

We used a container of 50 liters and a balloon of about 500 liters. When the gas flows out of the pressurized container, it slowly fills the balloon and when the pressure of the container reaches 1.5 bar, the compressor takes the gas from the balloon and fills the pressurized container.

The compressor worked in periods of about 2 minutes each 5 minutes, so that it remained within the specifications of the motor. The compressor used is a simple commercial membrane compressor. The scheme of the compressor is shown in Fig.3(b). We replaced the factory tubes by PVC pipings with pressure fittings. The compressor is activated by a pressure sensor installed in the container, that starts the motor at 1.5 bars and stops it at 4 bar. In addition, we installed two electrovalves that allow evacuating the overpressure in the exit of the compressor without emptying the container. This is needed as membrane compressors do not start with an overpressure at the exit. At the exit of the container, a flow regulator allows to establish the needed flow of gas through the circuit. A small filtering unit serves to remove impurities from the gas. A safety valve at 6 bar provides security for overpressure.

2.3 Sliding seal assembly and inner vacuum chamber for rapid sample exchange

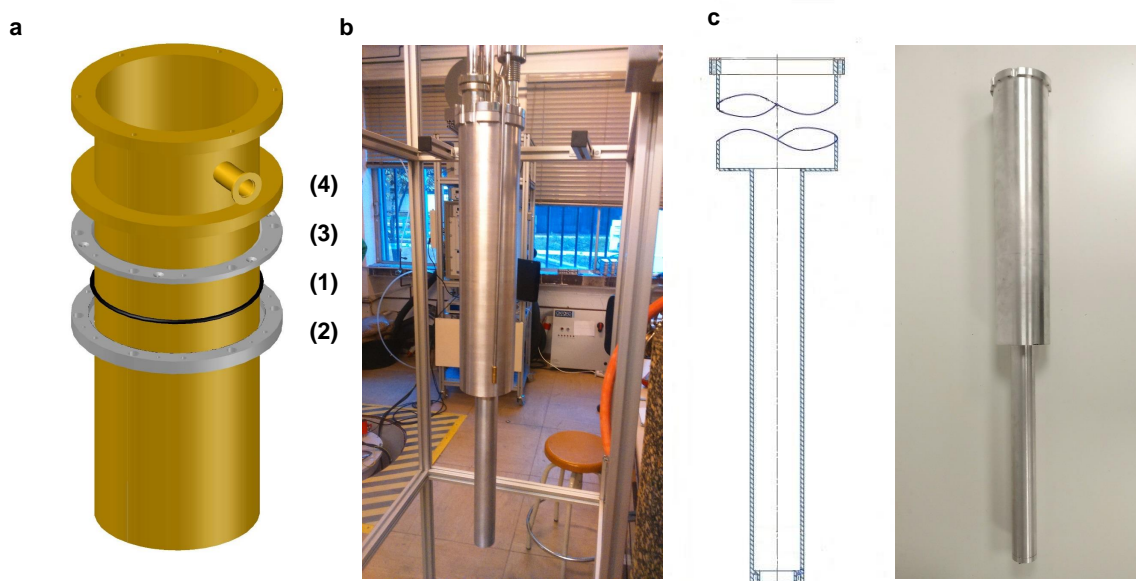


Figure 4: (a) Schematics of the sliding seal assembly consisting of an O-ring (1) and a bottom stainless steel piece (2) and a cap (3). The bottom part has threaded holes that allow fixing (3) in such a way as to compress the O-ring and additional holes allowing to fix the insert into the cryostat. (b) Photo of the insert outside and loaded into the cryostat. (c) Schematic drawing and photo of an Aluminium inner vacuum chamber.

Often, a refrigerator with diameter well above 100 mm has to be inserted into liquid Helium (see e.g. [19]). This can be made by a system that allows precooling the assembly using Helium gas. To this end, the cryostat needs a cylindrical enclosure with a Helium recovery port and a sliding seal that blocks the vertical motion of the cryostat, particularly if an overpressure appears in the cryostat, and is at the same time leak tight to avoid loss of Helium gas. A solution adopted by several companies is to use a radial shaft seal with a blocking device such as a spring[20]. The shaft seal is comprised by a plastic ring that is compressed using a spring. The shaft seal, however, is designed to be leak tight in presence of a fluid such as the lubricant of a shaft. When using it in a sliding seal assembly there is no lubricant and it can easily leak. To avoid leakage we have designed an O-Ring assembly comprised of a O-ring clamped using a stainless steel ring. The ring includes bolts to allow to fix the refrigerator to the cryostat and avoid motion if an overpressure is built up in the Dewar. It also includes an assembly that allows fixing the position of the insert when a overpressure is built up in the Dewar or when the Dewar is pumped.

Other solutions include the use of bellows and a seal assembly[21, 22]. These allow for lateral and rocking motions appearing during transport of a large cryostat assembly. In any event, these occupy considerably more space, hampering the usual operation when the cryostat is cold.

The bottom part of a dilution refrigerator contains the experiment and is located inside a vacuum can, often termed the inner vacuum chamber, as opposed to the 'outer' vacuum chamber of the containing Dewar. The inner vacuum chamber is most often made in stainless steel and

is partly immersed into liquid Helium. The vacuum chamber connects the upper flange of the chamber to the liquid Helium bath. The upper flange needs usually to be as cold as possible, to facilitate condensation of cryogens that arrive from room temperature through tubing. To this end, one solution is to weld a copper tube to the outside of the stainless steel inner vacuum chamber. This considerably increases the weight of the inner vacuum chamber and makes its cooling from room temperature costly in terms of evaporated liquid Helium.

To reduce the weight of the inner vacuum chamber, we have built a chamber out of Aluminum. The chamber comprises an indium seal assembly on the top to connect the chamber to the dilution insert. The chamber also comprises an indium seal on its bottom. The weight of the chamber is reduced by a factor of three, which correspondingly reduces the cryogenic liquid required for cooling or decreases the time needed to cool down again. The thermal conductivity of Aluminum alloy is about an order of magnitude higher than the thermal conductivity of stainless steel. This allows for effective cooling of the upper flange.

2.4 Acknowledgments

Work was supported through grant numbers FIS2017-84330-R, MDM-2014-0377 and RYC-2014-16626 of AEI, by the Comunidad de Madrid through program Nanofrontmag-CM (S2013/MIT-2850), by the European Research Council PNICTEYES grant agreement no. 679080 and by EU COST CA16218 Nanocohybri. Authors acknowledge the staff of segainvex at UAM.

3. Design files

3.1 Design Files Summary

Design file-name	File type	Open source license	Location of the file
Design file 1	Figure	CERN OHL	Available with the article, see Fig.1
DesignFile2 CAD and DesignFile2 STL	CAD plots	CERN OHL	https://osf.io/e6k7r/
DesignFile3 CAD and DesignFile3 STL	CAD plots	CERN OHL	https://osf.io/e6k7r/

Design file 1: Connecting schemes for the manual pump.

Design file 2: Design files, including all details, of the heat exchanger, of the scheme of the modifications of the compressor and of the parts made for the compressor.

Design file 3: Design files, including all details, of the sliding seal assembly and of the inner vacuum chamber.

4. Bill of materials

Designator	Component	Number	Cost per unit currency	Total cost	Source of materials	Material type
1	Manual pump	1	20\$	20\$	Usual sport or camping market	manual pump with 5 liter volume at each stroke and entry and exit ports
2	Rubber tubes and fittings	2	5\$	10\$	Tubes for the recovery	–

Designator	Component	Number	Cost per unit currency	Total cost	Source of materials	Material type
1	Compressor Stayer 70	1	280\$	280\$	Building market	High pressure machine
2	Solenoid valve SMC position normally open	2	50\$	100\$	Compressed air products	Valves
3	Presostat	1	50\$	50\$	Compressed air products	Sensor and valves
4	Balloon	1	1000\$	1000\$	Balloon made of quality rubber	Gas recipient
5	Fittings and pipes	–	70\$	70\$	–	–
Designator	Component	Number	Cost per unit currency	Total cost	Source of materials	Material type
1	Machining hardware	1	200\$	200\$	Hardware to machine the shown items	–

5. Build instructions

5.1 Manual pump

- Remove grease from the entry and exit ports of the manual pump.
- Use Araldite to make the pressure reading leak tight, putting epoxy on the juncture between the manometer and the tube.
- Use Araldite or other hard epoxy to glue connectors for flexible rubber tubes.
- Connect the entry of the manual pump to the recovery line and the exit to the transport Dewar.
- When pumping, take care not to build an overpressure in the transport Dewar. This is unlikely, as all transport Dewars have adequate safety valves, but it should nevertheless be noted.

5.2 Closed cycle circuit of cold Helium

- Substitute PVC tubes of the compressor with copper pipes, $\Phi = 6$ mm. Use conical copper fittings and make the required adaptations between tube diameters.
- Replace the pressostat of the compressor with a type Danfoss Kp35 060 pressostat or similar, a security valve (8 bar) and a filter with flow controller and valve.
- Replace the solenoid valves of the compressor with valves that are not leaking. Use SMC normally open valves. Be careful not to substitute the one-way valve at the compressor.
- Block the draining system of the tank with a screw.
- Seal with Araldite or similar all threads of the compressor.
- Connect the two entry ports of the compressor to a single line using a T.

5.3 Sliding seal assembly

- Make the hardware as shown in the design files.
- Screw the IVC using Indium seal as usual. Be careful and maintain a smooth aluminum surface to avoid leaks.

6. Operation instructions

- Usual precautions when using high pressure equipment. The manual pump does not require specific precautions. The compressor requires proper setting of security valves.
- When inserting the cryostat into liquid Helium always hold the cryostat by a thread so that it is not pushed out in presence of overpressure at the Helium recovery lines.

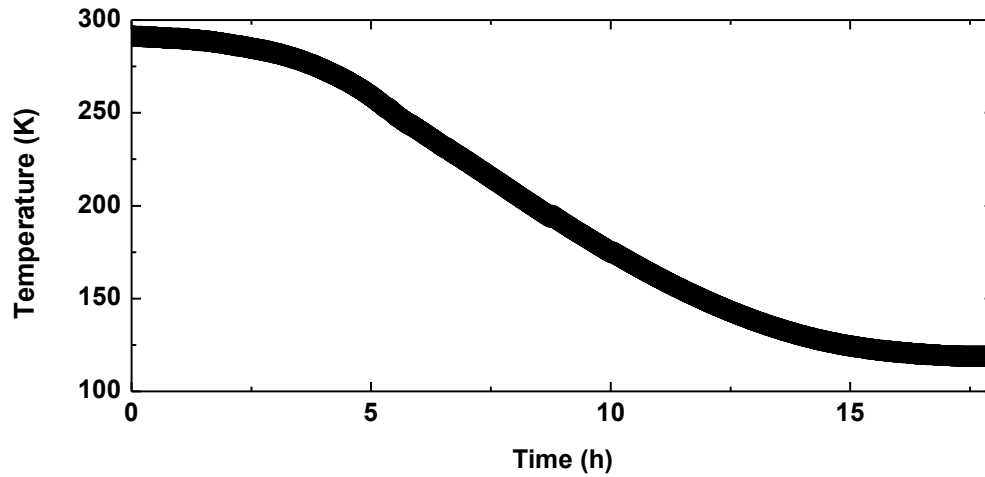


Figure 5: Temperature versus time taken during an overnight cool down of a 10 T superconducting magnet using the cold Helium closed-cycle system. During the morning, Helium was transferred with the same evaporation rate as found when using liquid nitrogen in the cryostat.

7. Validation and characterization

- We have observed that we can have transfer rates as high as 100 liters an hour operating a few times the manual pump. When the Helium transfer Dewar is nearly empty, the pump must be operated a few more times, as the volume needed to fill is much larger. We have used the manual pump in 50 liter and 100 liter transport Dewars. We do not expect major changes when using it in a 200 liter transport Dewar.
- In Fig.5 we show a typical cool down procedure to liquid nitrogen using the closed-cycle system. A bottle of about 60 liters of liquid nitrogen was emptied during the night.
- We have observed a decrease of the time needed to cool a cryostat by a factor of three when using the Aluminum inner vacuum chamber. In addition, the time needed to condense the mixture is reduced by a factor of two with respect to the use of stainless steel inner vacuum chambers.

8. Declaration of interest

Declarations of interest: none.

9. Human and animal rights

Does not apply.

References

- [1] Ray Radebaugh. *A Review of Pulse Tube Refrigeration*, pages 1191–1205. Springer US, Boston, MA, 1990.
- [2] Robert A. Ackermann, Kenneth G. Herd, and William E. Chen. *Advanced Cryocooler Cooling for MRI Systems*, pages 857–867. Springer US, Boston, MA, 2002.
- [3] Timothy John HughesKeith White. Cooling apparatus comprising a thermal interface and method for recondensing a cryogen gas, 2008. US20080216486A1.
- [4] S. H. Pan, E. W. Hudson, and J. C. Davis. The refrigerator based very low temperature scanning tunneling microscope. *Review of Scientific Instruments*, 70(2):1459–1463, 1999.
- [5] Young Jae Song, Alexander F. Otte, Vladimir Shvarts, Zuyu Zhao, Young Kuk, Steven R. Blankenship, Alan Band, Frank M. Hess, and Joseph A. Stroscio. Invited review article: A 10 mk scanning probe microscopy facility. *Review of Scientific Instruments*, 81(12):121101, 2010.
- [6] R. H. M. Smit, R. Grande, B. Lasanta, J. J. Riquelme, G. Rubio-Bollinger, and N. Agraït. A low temperature scanning tunneling microscope for electronic and force spectroscopy. *Review of Scientific Instruments*, 78(11):113705, 2007.
- [7] Marcus Liebmann, Jan Raphael Bindel, Mike Pezzotta, Stefan Becker, Florian Muckel, Tjorven Johnsen, Christian Saunus, Christian R. Ast, and Markus Morgenstern. An ultrahigh-vacuum cryostat for simultaneous scanning tunneling microscopy and magneto-transport measurements down to 400 mk. *Review of Scientific Instruments*, 88(12):123707, 2017.
- [8] S. Misra, B. B. Zhou, I. K. Drozdov, J. Seo, L. Urban, A. Gyenis, S. C. J. Kingsley, H. Jones, and A. Yazdani. Design and performance of an ultra-high vacuum scanning tunneling microscope operating at dilution refrigerator temperatures and high magnetic fields. *Review of Scientific Instruments*, 84(10):103903, 2013.
- [9] Maximilian Assig, Markus Etzkorn, Axel Enders, Wolfgang Stiepany, Christian R. Ast, and Klaus Kern. A 10 mk scanning tunneling microscope operating in ultra high vacuum and high magnetic fields. *Review of Scientific Instruments*, 84(3):033903, 2013.
- [10] T. Quaglio, F. Dahlem, S. Martin, A. Gérardin, C. B. Winkelmann, and H. Courtois. A subkelvin scanning probe microscope for the electronic spectroscopy of an individual nano-device. *Review of Scientific Instruments*, 83(12):123702, 2012.
- [11] T. Machida, Y. Kohsaka, and T. Hanaguri. A scanning tunneling microscope for spectroscopic imaging below 90 mk in magnetic fields up to 17.5 t. *Review of Scientific Instruments*, 89(9):093707, 2018.
- [12] Nathaniel Craig and Ted Lester. Hitchhiker’s guide to the dilution refrigerator, 2004.
- [13] N. Howard Balshaw. Practical cryogenics: An introduction to laboratory cryogenics. ISBN-10 0952759403.
- [14] H. Suderow, I. Guillamon, and S. Vieira. Compact very low temperature scanning tunneling microscope with mechanically driven horizontal linear positioning stage. *Review of Scientific Instruments*, 82(3):033711, 2011.

- [15] J. A. Galvis, E. Herrera, I. Guillamón, J. Azpeitia, R. F. Luccas, C. Munuera, M. Cuenca, J. A. Higuera, N. Díaz, M. Pazos, M. García-Hernandez, A. Buendía, S. Vieira, and H. Suderow. Three axis vector magnet set-up for cryogenic scanning probe microscopy. *Review of Scientific Instruments*, 86(1):013706, 2015.
- [16] Beat MrazChristoph NotzliErnst UngrichtDaniel Marek. Method and a device for precooling the helium tank of a cryostat, 1989. US5187938A.
- [17] L. Skrbek. Turbulence in cryogenic helium. *Physica C: Superconductivity*, 404(1):354 – 362, 2004. Proceedings of the Third European Conference on Vortex Matter in Superconductors at Extreme Scales and Conditions.
- [18] J. R. Roebuck and H. Osterberg. The joule-thomson effect in helium. *Phys. Rev.*, 43:60–69, Jan 1933.
- [19] Peter Mile Cetnik Paul Kelly, David Massingham and Christopher David Wilkinson. Dilution refrigerator assembly, 2001. US Patent 6758059B2.
- [20] Ferdinand Schweitzer. Radial shaft seal with retaining device for a tension spring, 2010. US Patent 8668390B2.
- [21] Evangelos T. Laskaris. Support for cryostat penetration tube, 1984. US Patent 4526015A.
- [22] Evangelos T. Laskaris. Horizontal cryostat penetration insert and assembly, 1984. US Patent 4522034A.