

Next Generation Radiogenomics Sequencing for Prediction of EGFR and KRAS Mutation Status in NSCLC Patients Using Multimodal Imaging and Machine Learning Approaches

Isaac Shiri^{1,2}, Hassan Maleki^{1,3}, Ghasem Hajianfar^{1,2}, Hamid Abdollahi^{1,4}, Saeed Ashrafinia^{5,6}, Mathieu Hatt⁷, Mehrdad Oveisi^{1,8}, Arman Rahmim^{5,9,10}

1. Department of Biomedical and Health Informatics, Rajaie Cardiovascular Medical and Research Center, Iran University of Medical Science, Tehran, Iran
2. Research Center for Molecular and Cellular Imaging, Tehran University of Medical Sciences, Tehran, Iran
3. Department of Computer Science and Engineering, Shahid Beheshti University, Tehran, Iran
4. Department of Medical Physics, School of Medicine, Iran University of Medical Science, Tehran, Iran
5. Department of Radiology and Radiological Science, Johns Hopkins University, Baltimore MD, USA
6. Department of Electrical and Computer Engineering, Johns Hopkins University, Baltimore MD, USA
7. INSERM, UMR 1101, LaTIM, Univ Brest, F-29238, Brest, France
8. Department of Computer Science, University of British Columbia, Vancouver BC, Canada
9. Departments of Radiology and Physics & Astronomy, University of British Columbia, Vancouver BC, Canada
10. Department of Integrative Oncology, BC Cancer Research Centre, Vancouver BC, Canada

First Author: Isaac Shiri, MSc

Department of Biomedical and Health Informatics, Rajaie Cardiovascular Medical and Research Center, Iran University of Medical Science, Tehran, Iran

Phone: +989190245882

Email: Isaac.sh92@gmail.com

Corresponding Author: Arman Rahmim, PhD

Associate Professor of Radiology and Physics, University of British Columbia
Senior Scientist & Provincial Medical Imaging Physicist, BC Cancer Agency
BC Cancer Research Center

675 West 10th Ave

Office 5-114

Vancouver, BC, V5Z 1L3

<http://rahmimlab.com>

Phone: 604-675-8262

Email: arman.rahmim@ubc.ca

Conflicts of Interest: The authors declare no potential conflicts of interest.

Abstract

Purpose

Considerable progress has been made in assessment and management of NSCLC patients based on mutation status in the EGFR and KRAS. At the same time, NSCLC management through KRAS and EGFR mutation profiling faces some challenges. In the present work, we aimed to evaluate a comprehensive radiomics framework that enabled prediction of EGFR and KRAS mutation status in NSCLC cancer patients based on low dose CT, diagnostic CT, PET modalities radiomic features and machine learning (ML) algorithms.

Methods

Our study involved NSCLC cancer patient with 186 PET, and 175 low dose CT and CTD images. More than twenty thousand radiomic features from different image-feature sets were extracted. Conventional clinically used standard uptake value (SUV) parameters were also obtained for PET images. Feature value was normalized to obtain Z-scores, followed by student t-test students for comparison, high correlated features were eliminated and the False discovery rate (FDR) correction were performed and q-value were reported for univariate analysis. Six feature selection methods and twelve classifiers were used to predict gene status in patient. We performed 10-fold cross-validation for model tuning to improve robustness and all model evaluation was reported on independent validation sets (68 patients). The mean area under the receiver operator characteristic (ROC) curve (AUC) was obtained for performance evaluation.

Results

The best predictive power of conventional PET parameters was achieved by SUVpeak (AUC: 0.69, P-value = 0.0002) and MTV (AUC: 0.55, P-value = 0.0011) for EGFR and KRAS, respectively. Univariate analysis of extracted radiomics features improved prediction power up to AUC: 0.75 (q-value: 0.003, Short Run Emphasis feature of GLRLM from LOG preprocessed image of PET with sigma value 1.5) and AUC: 0.71 (q-value 0.00005, The Large Dependence Low Gray Level Emphasis from GLDM in LOG preprocessed image of CTD sigma value 5) for EGFR and KRAS, respectively. Furthermore, the machine learning algorithm improved the prediction power up to AUC: 0.82 for EGFR (LOG preprocessed of PET image set with sigma 3 with VT feature selector and SGD classifier) and AUC: 0.83 for KRAS (CT image set with sigma 3.5 with SM feature selector and SGD classifier).

Conclusion

Our findings demonstrated that non-invasive and reliable radiomics analysis can be successfully used to predict EGFR and KRAS mutation status in NSCLC patients. We demonstrated that radiomic features extracted from different image-feature sets could be used for EGFR and KRAS mutation status prediction in NSCLC patients, and showed that they have more predictive power than conventional imaging parameters.

Key Words: Radiogenomics, PET/CT, Machine Learning, NSCLC, KRAS, EGFR

Introduction

Considerable progress has been made recently in the assessment and management of non-small cell lung cancer (NSCLC) based on mutation status in the epidermal growth factor receptor (EGFR) and Kirsten rat sarcoma viral oncogene (KRAS) genes (1). Ongoing studies on molecular cancer profiling have revealed that EGFR and KRAS are involved in the occurrence, development, invasion, and metastasis of NSCLC (2). Moreover, studies have indicated that mutations in KRAS and EGFR are considered as first lines for clinical decision making in NSCLC treatment and outcome improvement (3). Furthermore, recent studies have identified that NSCLC patients with mutant KRAS tumors fail to benefit from adjuvant chemotherapy, and their disease does not respond to EGFR inhibitors (4). In addition, it was found that outcomes and patterns of failure in genotypic subgroups of NSCLC patients, based on mutations in EGFR or KRAS, can inform the design of future trials integrating targeted therapies (5).

Although KRAS and EGFR mutation profiling is critical in NSCLC management, some studies have raised issues with this approach (6, 7). First, such mutation status captures only a small degree of tumor heterogeneity and does not provide a complete picture for the assessment of tumor characteristics. Secondly, this method depicts low repeatability and is not feasible for all cases (8). Furthermore, the method suffers from invasiveness and patients discomfort (9).

As a recently developed paradigm of advanced medical image quantification, radiomics has garnered significant interests given its cost-effectiveness and reliability to characterize tumor heterogeneity, and has enabled improved assessment of therapy response and prediction of molecular pathways (10-14). Accumulating evidence has identified several radiomic features extracted from computed tomography (CT), magnetic resonance (MR) or positron emission tomography (PET) images as highly correlated with genomic parameters in several cancers (8, 15-18). For NSCLC patients, radiomics studies have shown several CT image-features can predict mutation status in EGFR and KRAS (8). For example, features such as size, edge, lucency and homogeneity extracted from CT images could identify EGFR mutation status (19). Also, some diffusion-weighted (DW) MR image-intensity histogram features, including mean, skewness, and 10th and 90th percentiles, have been shown to predict EGFR mutation in lung adenocarcinoma (20). On the other hand, Velazquez *et al.* developed radiomics models based on CT image-features and

clinical parameters to distinguish between EGFR⁻ and EGFR⁺, and KRAS⁺ and KRAS⁻ (21). Liu *et al.* also evaluated the ability of CT image-features to predict EGFR mutation status in 298 surgically-resected peripheral lung adenocarcinomas in an Asian cohort of patients, and built a high performance predictive model by using multiple logistic regression algorithms (19). Zhang *et al.* also developed a radiogenomic model based on CT image-features to predict EGFR and KRAS mutations in lung adenocarcinoma patients (22).

Incorporating radiomics for a prediction study requires a multi-step process that involves reliable statistical analyses, such as feature selection and classification, to reduce over-fitting and to build robust predictive or prognostic models (23). A number of machine learning (ML) algorithms can provide robust means to identify a subset of features to combine into a multi-parametric model (24). Although several ML algorithms, alone or in combination, have been used in radiomics analysis for feature selection and classification, there is no “one fits all” approach as performance of various ML workflows have been shown to depend on application and/or type of data (25-27) Previous studies have tested cross-combination of different ML approaches, and have suggested distinct ML algorithms that depict high performance for feature selection and classification (24, 25, 28).

In the present study, we aimed to evaluate a comprehensive multimodal (Diagnostic CT, low dose CT and PET modality) radiomics, univariate analysis and machine learning framework to predict EGFR and KRAS mutation status in NSCLC cancer patients.

Material and methods

Radiomics analysis

Our radiomics analysis included the following seven steps (Fig. 1).

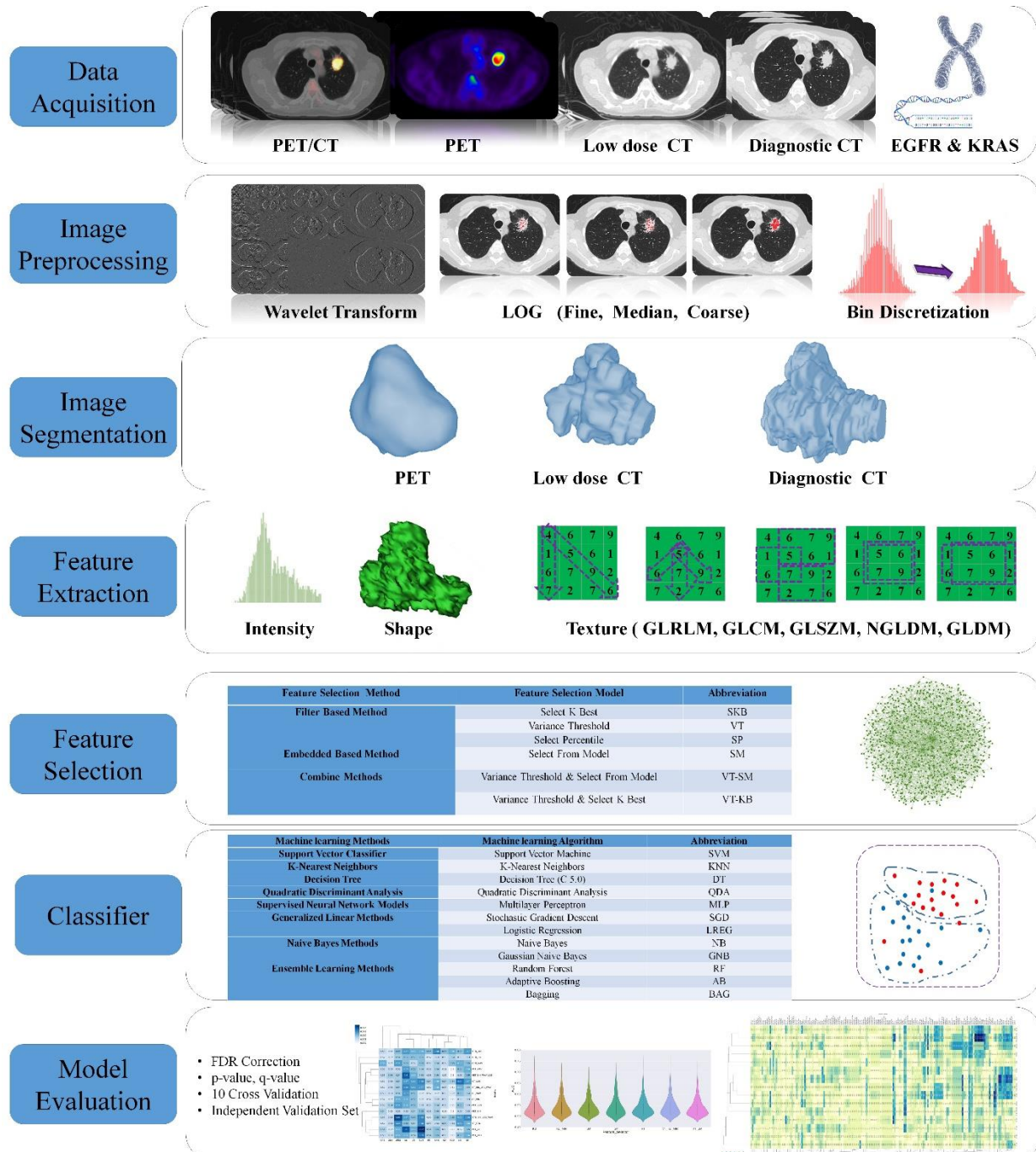


Figure 1. The radiogenomics framework employed in the present study

1- Data collection

Clinical characteristics of all patients were presented in Table 1.

Table 1. Clinical characterization of patients

Characteristic		PET	CT	CTD	Validation PET, CT, CTD
Patient NO.		186	175	175	68
Age, mean(SD)		68.6 (9.11)	69.2 (8.85)	69.2 (8.85)	69.1 (8.8)
Weight (Kg), mean(SD)		77.0 (18.5)	77.8 (18.6)	77.8 (18.6)	77.6 (18.6)
Sex (%)	Male	66.7	66.3	66.3	67.1
	Female	33.3	33.7	33.7	32.9
Smoking status (%)	Nonsmoker	21	21.7	21.7	21.5
	Former	61.3	61.7	61.7	61.0
	Current	17.7	16.5	16.5	17.3
Histology (%)	Adenocarcinoma	79.6	82.3	82.3	81.4
	Squamous cell carcinoma	18.3	16	16	16.8
	NOS (not otherwise specified)	2.2	1.7	1.7	1.8
T stage (%)	T1	36.1	38.2	38.2	37.1
	T2	30.6	28	28	29.3
	T3	11.3	9.1	9.1	9.6
	T4	3.8	2.9	2.9	3
N stage (%)	N0	66.7	63.4	63.4	65.3
	N1	8.1	6.9	6.9	7.2
	N2	9.7	9.7	9.7	10.2
M stage (%)	M0	82.3	78.3	78.3	80.8
	M1	2.2	1.7	1.7	1.8
AJCC Staging (%)	I	50	50.3	50.3	51.5
	II	19.4	16	16	16.8
	III	12.9	12	12	12.6
	III	2.2	1.7	1.7	1.8
Histo-pathological Grade (%)	G1	15.6	14.9	14.9	15
	G2	39.8	38.9	38.9	40.1
	G3	17.7	15.4	15.4	16.2
	other	11.3	10.8	10.8	11.4
Lymph vascular invasion (%)	Absent	68.3	67.4	67.4	69.5
	Present	11.3	10.3	10.3	10.8
EGFR mutation status (%)	Mutant	18.8	21.1	21.1	19.8
	Wild-type	61.8	62.9	62.9	63.5
	Not Available	18.8	16	16	16.8
KRAS mutation status (%)	Mutant	18.3	18.3	18.3	18.6
	Wild-type	61.3	65.1	65.1	64.1
	Not Available	20.4	16.6	16.6	17.4
PET SUV	SUVmax	7.5 (6.7)			7.1 (6.3)
	SULmax	6.1(5.2)			5.6 (4.9)
	SUVpeak	9.6 (6.5)			8.9 (6.3)
	SULpeak	7.4 (4.9)			6.8 (4.8)
Volume mean(SD)		29915 (48293)			19225 (40185)
Recurrence (%)		17.7	17.7	17.7	17.4
Time to Last Follow-up, mean(SD)		26.9 (21.7)	27.0 (22.1)	27.0 (22.1)	26.8 (22.2)
Death (%)		21	20	20	20.4

1-1- Imaging

This study was conducted on 211 NSCLC cancer patients with available imaging and genomic data. Imaging modalities, including diagnostic CT (CTD) and PET/CT (i.e., low-dose CT (CT) used for PET attenuation correction and the PET image) for all patients were obtained (29-32). All images were acquired with two vendor of imaging system including; GE and SIEMENS. For CTD, 73% of data were acquired with GE Discovery CT750HD and 27% with SIEMENS. Slice Thickness for diagnostics were 0.625–3 mm (median: 1.5 mm) and an X-ray tube current of 124–699 mA (mean 220 mA) at 80–140 kVp (mean 120 kVp), and Spiral Pitch Factor 0.9-1. A GE Discovery D690 PET/CT and GE Discovery PET/CT scanner were used for PET/CT scanning. 18F-FDG Dose and uptake time were 138.90–572.25 MBq (mean: 309.26 MBq) and 23.08–128.90 min (mean: 66.58 minutes), respectively. Each bed position was 1–5-minute acquisition, and iterative Ordered Subset Expectation Maximization (OSEM) reconstruction was used for PET image reconstruction. Slice Thickness for low dose CT were set as 3-5mm, mA Range 50-400, KVP 120-140 and Spiral Pitch Factor 0.9-1. Due to certain problems including image noise, artifacts, absence of images for some sequences and image mis-segmentation, some patients were excluded, and finally 186 PET, as well as 175 low dose CT and CTD images were segmented.

1-2- Genomics

Tumor samples were excised and cutted by surgeon with slice thickness of 3-5 mm and were frozen for 30 minutes. For single nucleotide mutation detection, SNaPshot technology based on dideoxy single-base extension of oligonucleotide primers was performed after multiplex polymerase chain reaction (PCR). To detect EGFR mutations, Exons including 18, 19, 20 and 21 were assessed. In addition, Exon 2 Positions 12 and 13 were examined with amino acid substitution for missense KRAS mutations.

2- Image pre-processing

Prior to feature extraction, all images were pre-processed with interpolation to isotropic voxel spacing to be rotationally invariant, for feature extraction to allow comparison between image data from different samples and scanner and center. Re-sampling to $1 \times 1 \times 1$, $1 \times 1 \times 3$ and $3 \times 3 \times 3$ were performed for CTD, CT and PET images respectively. Laplacian of Gaussian (LOG), wavelet decomposition (WAV) and discretized into 64 bins (BIN64) preprocessing were performed to generate different set of features. For LOG filter, different sigma values were used to extract fine, medium and coarse features; specifically, they ranged from 0.5 to 5 with 0.5 steps. Wavelet filtering, yields 8 decompositions per level (all possible combinations of applying either a High or a Low pass filter in each of the three dimensions including HHH, HHL, HLH, HLL, LHH, LHL, LLH and LLL). Preprocessing steps (including discretization, LOG and wavelet) were performed on all intensity, histogram and textural features.

3- Image segmentation

All segmentations in PET images were performed using OSIRIX[®] (33). Lesions were delineated manually on PET imaged. CT images were segmented via automatic region growing using 3D-Slicer (34), and edited and verified by an experienced radiologist. In total we provide the 186 PET, and 175 low dose CT and CTD segmentation region.

4- Feature extraction

In the next step, several features from different feature-classes were extracted. The classes include first order statistics (19 FOS features), shape-based (16 Shape features), gray level co-occurrence matrix (GLCM 23 features), gray level run length matrix (16 GLRLM features), gray level size zone matrix (16 GLSZM features), neighboring gray tone difference matrix (5 NGTDM features), and gray level dependence matrix (14 GLDM features). Image feature extraction was performed by the open-source python library PyRadiomics (35). Radiomics features calculated by this package are in compliance with feature definitions as described by the image biomarker standardization initiative (IBSI) that ensures the harmonization and reproducibility of calculated radiomic features, and thus facilitates the reproducibility of this study (35, 36). More detail on the feature classes is provided in Table 2. In addition to the radiomic features, we also extracted conventional clinical PET biomarkers including metabolic tumor volume (MTV), standard uptake values (SUV_{max} , SUV_{peak}), including normalization to lean body mass (SUL_{max} , SUL_{peak}).

Table 2. Radiomic features extracted from different modality

First Order Statistics (FOS)	Gray Level Co-occurrence Matrix (GLCM)	Gray Level Run Length Matrix (GLRLM)
1. Energy	1. Autocorrelation	1. Short Run Emphasis (SRE)
2. Total Energy	2. Joint Average	2. Long Run Emphasis (LRE)
3. Entropy	3. Cluster Prominence	3. Gray Level Non-Uniformity (GLN)
4. Minimum	4. Cluster Shade	4. Gray Level Non-Uniformity Normalized (GLNN)
5. 10th percentile	5. Cluster Tendency	5. Run Length Non-Uniformity (RLN)
6. 90th percentile	6. Contrast	6. Run Length Non-Uniformity Normalized (RLNN)
7. Maximum	7. Correlation	7. Run Percentage (RP)
8. Mean	8. Difference Average	8. Gray Level Variance (GLV)
9. Median	9. Difference Entropy	9. Run Variance (RV)
10. Interquartile Range	10. Difference Variance	10. Run Entropy (RE)
11. Range	11. Joint Energy	11. Low Gray Level Run Emphasis (LGLRE)
12. Mean Absolute Deviation (MAD)	12. Joint Entropy	12. High Gray Level Run Emphasis (HGLRE)
13. Robust Mean Absolute Deviation (rMAD)	13. Informal Measure of Correlation (IMC) 1	13. Short Run Low Gray Level Emphasis (SRLGLE)
14. Root Mean Squared (RMS)	14. Informal Measure of Correlation (IMC) 2	14. Short Run High Gray Level Emphasis (SRHGLE)
15. Standard Deviation	15. Inverse Difference Moment (IDM)	15. Long Run Low Gray Level Emphasis (LRLGLE)
16. Skewness	16. Inverse Difference Moment Normalized (IDMN)	16. Long Run High Gray Level Emphasis (LRHGLE)
17. Kurtosis	17. Inverse Difference (ID)	
18. Variance	18. Inverse Difference Normalized (IDN)	Gray Level Dependence Matrix (GLDM)
19. Uniformity	19. Inverse Variance	1. Small Dependence Emphasis (SDE)
	20. Maximum Probability	2. Large Dependence Emphasis (LDE)
	21. Sum Average	3. Gray Level Non-Uniformity (GLN)
	22. Sum Entropy	4. Dependence Non-Uniformity (DN)
	23. Sum of Squares	5. Dependence Non-Uniformity Normalized (DNN)
		6. Gray Level Variance (GLV)
		7. Dependence Variance (DV)
		8. Dependence Entropy (DE)
		9. Low Gray Level Emphasis (LGLE)
		10. High Gray Level Emphasis (HGLE)
		11. Small Dependence Low Gray Level Emphasis (SDLGLE)
		12. Small Dependence High Gray Level Emphasis (SDHGLE)
		13. Large Dependence Low Gray Level Emphasis (LDLGLE)
		14. Large Dependence High Gray Level Emphasis (LDHGLE)
		Neighboring Gray Tone Difference Matrix (NGTDM)
		1-Coarseness
		2-Contrast
		3-Busyness
		4-Complexity
		5- Strength
Shape and Morphological Features	Gray Level Size Zone Matrix (GLSZM)	
1. Volume	1. Small Area Emphasis (SAE)	
2. Surface Area	2. Large Area Emphasis (LAE)	
3. Surface Area to Volume ratio	3. Gray Level Non-Uniformity (GLN)	
4. Sphericity	4. Gray Level Non-Uniformity Normalized (GLNN)	
5. Compactness 1	5. Size-Zone Non-Uniformity (SZN)	
6. Compactness 2	6. Size-Zone Non-Uniformity Normalized (SZNN)	
7. Spherical Disproportion	7. Zone Percentage (ZP)	
8. Maximum 3D diameter	8. Gray Level Variance (GLV)	
9. Maximum 2D diameter (Slice)	9. Zone Variance (ZV)	
10. Maximum 2D diameter (Column)	10. Zone Entropy (ZE)	
11. Maximum 2D diameter (Row)	11. Low Gray Level Zone Emphasis (LGLZE)	
12. Major Axis	12. High Gray Level Zone Emphasis (HGLZE)	
13. Minor Axis	13. Small Area Low Gray Level Emphasis (SALGLE)	
14. Least Axis	14. Small Area High Gray Level Emphasis (SAHGLE)	
15. Elongation	15. Large Area Low Gray Level Emphasis (LALGLE)	
16. Flatness	16. Large Area High Gray Level Emphasis (LAHGLE)	

5- Univariate analysis

For univariate analysis, each feature value was normalized to obtain Z-scores, followed by student t-test students for comparison. A p-value of <0.05 was used as a criterion for statistically significant results. Spearman correlation between features were performed to eliminate high correlated features and the False discovery rate (FDR) correction were assessed and q-value were reported. Statistical analysis was performed in R 3.5.1 (using ‘pROC’ and ‘stats’ packages). We also studied the prediction performance of conventional clinical PET parameters (SUV_{max} , SUV_{peak} , SUL_{max} , SUL_{peak} and MTV).

6- Feature selection

We implemented 6 different feature selections methods in our framework and compared their performances (see Table 3). We implement 6 feature selection including Filter Based Method (Select K Best (SKB), Variance Threshold (VT) and Select Percentile (SP)), Embedded Based Method (Select from Model (SM)) and Combine Methods (Variance Threshold & Select from Model (VT-SM) and Variance Threshold & Select K Best (VT-KB)).

Table 3. Feature selector algorithm characteristics

Feature Selection Method	Feature Selection Model	Abbreviation	Algorithm of action
Filter Based Method	Select K Best	SKB	This method scores the features using a chi square function and then "removes all but the k highest scoring features.
	Variance Threshold	VT	It removes all features whose variance doesn't meet some threshold. By default, it removes all zero-variance features, i.e. features that have the same value in all sample. In this method, we used 0.8 threshold.
	Select Percentile	SP	It removes all but a user-specified highest scoring percentage of features using common univariate statistical tests for each feature: false positive rate or family wise error
Embedded Based Method	Select From Model	SM	It is a meta-transformer that can be used along with any estimator that has a coefficient or feature importance's attribute after fitting. The features are considered unimportant and removed, if the corresponding coefficient or feature importance's values are below the provided threshold parameter. Apart from specifying the threshold numerically, there are built-in heuristics for finding a threshold using a string argument. Available heuristics are "mean", "median" and float multiples of these like "0.1*mean".
Combine Methods	Variance Threshold & Select From Model	VT-SM	Combine of Variance Threshold & Select From Model
	Variance Threshold & Select K Best	VT-KB	Combine of Variance Threshold & Select K Best

7- Classifier

We implemented and compared 12 classifiers (Table 4, and detail of each classifier are provided in supplemental Table 1). Different method of classifications including generalized linear models (logistic regression and stochastic gradient descent), naive Bayes models (naive Bayes and Gaussian naive Bayes), nearest neighbor's model (k-nearest neighbors), decision trees model (C5.0), quadratic discriminant analysis model (QDA), support vector machines (SVC), supervised neural network models (multi-layer perceptron) and ensemble learning methods (adaptive boost, bagging and random forest) were used in this study.

Table 4. Classifier methods and abbreviation

Machine learning Methods	Machine learning Algorithm	Abbreviation	Algorithm of action
Support Vector Classifier	Support Vector Machine	SVM	Data items are plotted as a point in n-dimensional space (n is number of features) with the value of each feature being the value of a particular coordinate and classification is performed by finding the hyper-plane that differentiate the two classes.
K-Nearest Neighbors	K-Nearest Neighbors	KNN	As a non-parametric, lazy learning algorithm, KNN algorithm is based on feature similarity and an object is classified by a majority vote of its neighbors, with the object being assigned to the class most common among its k nearest neighbors.
Decision Tree	Decision Tree (C 5.0)	DT	As a decision support tool, uses a tree-like graph of decisions and their possible consequences, including chance event outcomes, resource costs, and utility, builds a model of decisions.
Quadratic Discriminant Analysis	Quadratic Discriminant Analysis	QDA	QDA separate measurements of two classes of objects with a quadratic decision boundary, generated by fitting class conditional densities to the data. It is a more general version of the linear classifier.
Supervised Neural Network Models	Multilayer Perceptron	MLP	MLP is a class of feed forward artificial neural network which optimizes the log-loss function using stochastic gradient descent (SGD).
Generalized Linear Methods	Stochastic Gradient Descent	SGD	Works based on the differentiable objective function optimization which is a stochastic approximation of gradient descent optimization. This classifier implements a plain SGD learning routine which supports different loss functions and penalties for classification.
	Logistic Regression	LREG	The probabilities describing the possible outcomes of a single trial are modeled using a logistic function.
Naive Bayes Methods	Naive Bayes	NB	As a probabilistic classifier, works based on applying Bayes' theorem with strong (naive) independence assumptions between the features.
	Gaussian Naive Bayes	GNB	It is a generalization of Naive Bayes Networks, which are a special case of probabilistic networks that allows treating continuous variables. This method computes conditional class probabilities and then predicts the most probable class of a vector of training data X, according to sample data D.
Ensemble Learning Methods	Random Forest	RF	It is an ensemble learning approach which is essentially a collection of decision trees, where each tree is slightly different from the others. This algorithm reduces the amount of overfitting by averaging the result of all decision trees.
	Adaptive Boosting	AB	This ensemble learning method, consists of very simple base classifiers (weak learners). This approach focuses on the training samples that are hard to classify, that is, to let the weak learners subsequently learn from misclassified training samples to improve the performance of the ensemble.
	Bagging	BAG	BAG fits base classifiers each on random subsets of the original dataset and then aggregate their individual predictions (by voting) to form a final prediction. The base estimator to fit on random subsets of the dataset is a decision tree model.

8- Model evaluation

All our analyses, including feature selection and classification, were performed using an in-house developed python framework in open-source python library Scikit-Learn (37). Cross-validation (CV) was applied to models and tuning of models were performed on testing set of 10-fold CV and this process were performed 20 times to get highest stable results. Furthermore, models evaluation was performed on independent validation sets (65 patients). The predictive power of all features was investigated using the area under the receiver operator characteristic (ROC) curve (AUC). A heatmap was prepared to compare different developed models. Mean \pm standard deviation (SD) of AUC of all classifier and feature selection algorithms were calculated and depicted. The cross combination of feature selection and classification methods were depicted as a heatmap (mean AUC value in cross validation).

Results

Univariate analysis

- **EGFR**

Results of EGFR mutation status prediction after elimination of high correlated feature in each image sets and FDR correction were presented in supplemental Table 2 (only q-value <0.05). About ninety features from all image data set had significant correlation with EGFR mutation status (q-value < 0.05) with mean AUC 0.67 ± 0.05 (0.52-0.75). Highest prediction performance was achieved by Short Run Emphasis feature of GLRLM (AUC: 75, q-value: 0.003) and Run Variance feature of GLRLM (AUC: 75, q-value: 0.016) from LOG preprocessed image of PET with sigma value 1.5 and 1, respectively. Also, Small Area Emphasis feature and Size Zone Non Uniformity Normalized from GLSZM of CTD data set which preprocessed by LLL of wavelet and LOG with sigma 2 respectively (AUC: 75, q-value: 0.001). More detail of univariate analysis in EGFR mutation status prediction including AUC, p-value and q-value were presented in supplementary heatmap 1-3, respectively.

- **KRAS**

Results of KRAS mutation status prediction after elimination of high correlated feature in each image sets and FDR correction were presented in supplemental Table 3 (only q-value <0.05). Fourteen features from all image data set were significantly correlated with KRAS mutation status (q-value < 0.05) with mean AUC 0.67 ± 0.02 (0.53-0.71). The Large Dependence Low Gray Level Emphasis from GLDM in LOG (sigma: 5) preprocessed image of CTD with AUC 0.71 (q-value: 0.00005) had highest performances. Also, the Minimum feature from first order of LOG preprocessed in sigma 4 (AUC: 71, q-value: 0.013), 4.5 (AUC: 71, q-value: 0.01) and 5 (AUC: 70, q-value: 0.014) showed high performances in KRAS mutation status prediction. More detail of univariate analysis in KRAS mutation status prediction including AUC, p-value and q-value were presented in supplementary heatmap 4-6, respectively.


Supplemental Table 2. EGFR mutation status perdition after elimination of high correlated feature in each image sets and FDR correction were presented as AUC, p-value and q-value

Image_Sets	Preprocessed	Filter	Type	Features	AUC	p.value	q-value			
PET	BIN 64	BIN64	FO	10Percentile	0.64	0.003	0.024			
			GLDM	DV	0.71	0.007	0.027			
	LOG	Sigma: 0.5	firstorder	Mean	0.68	0.000	0.001			
				90Percentile	0.63	0.001	0.003			
			glcm	Idn	0.67	0.000	0.002			
				MaximumProbability	0.73	0.001	0.002			
				Imc1	0.63	0.003	0.008			
			gldm	DependenceVariance	0.73	0.000	0.002			
				SmallDependenceLowGrayLevelEmphasis	0.54	0.017	0.037			
			glszm	SmallAreaEmphasis	0.70	0.000	0.002			
				GrayLevelNonUniformity	0.63	0.001	0.004			
				SizeZoneNonUniformity	0.71	0.002	0.006			
			Sigma: 1.0	firstorder	90Percentile	0.63	0.001	0.009		
					glcm	Imc2	0.70	0.004	0.016	
		glrlm		RunVariance	0.75	0.005	0.016			
		Sigma: 1.5	firstorder	Mean	0.68	0.000	0.000			
				90Percentile	0.62	0.001	0.003			
				Median	0.58	0.018	0.044			
			glrlm	ShortRunEmphasis	0.75	0.001	0.003			
		Sigma: 2.0	firstorder	Mean	0.68	0.000	0.001			
				90Percentile	0.60	0.002	0.006			
				Median	0.61	0.003	0.008			
				Mean	0.69	0.000	0.000			
			glszm	ZonePercentage	0.73	0.000	0.001			
				gldm	DependenceNonUniformityNormalized	0.73	0.000	0.000		
		Sigma: 3.0	firstorder	Mean	0.69	0.000	0.001			
				Median	0.64	0.000	0.001			
			glrlm	ShortRunEmphasis	0.73	0.000	0.001			
		Sigma: 3.5	firstorder	Median	0.67	0.000	0.000			
				Maximum	0.54	0.005	0.015			
			gldm	LargeDependenceEmphasis	0.73	0.001	0.006			
		Sigma: 4.5	firstorder	Median	0.69	0.000	0.000			
			gldm	LargeDependenceEmphasis	0.71	0.002	0.008			
		Sigma: 4	firstorder	Median	0.68	0.000	0.000			
				Maximum	0.52	0.016	0.049			
	glrlm		RunLengthNonUniformityNormalized	0.70	0.002	0.009				
	Sigma: 5.0	firstorder	Median	0.70	0.000	0.000				
		gldm	LargeDependenceEmphasis	0.69	0.003	0.017				
	WAV	WAV	HHH	gldm	DependenceNonUniformityNormalized	0.71	0.000	0.000		
				glszm	SmallAreaEmphasis	0.65	0.002	0.017		
			HHL	firstorder	Mean	0.61	0.010	0.039		
			HHL	gldm	DependenceNonUniformityNormalized	0.73	0.000	0.001		
			HHL	ngtdm	Strength	0.73	0.005	0.028		
			HLH	firstorder	Median	0.63	0.009	0.034		
				glrlm	RunLengthNonUniformityNormalized	0.73	0.000	0.003		
				ngtdm	Strength	0.71	0.008	0.034		
			HLL	glszm	SizeZoneNonUniformityNormalized	0.72	0.000	0.002		
			LHH	glrlm	RunPercentage	0.73	0.001	0.005		
			LHL	firstorder	Median	0.59	0.011	0.026		
				gldm	SmallDependenceEmphasis	0.73	0.000	0.002		
				glszm	LargeAreaHighGrayLevelEmphasis	0.66	0.002	0.010		
				ngtdm	Busyness	0.65	0.009	0.026		
			LLH	firstorder	Mean	0.66	0.001	0.005		
				glcm	InverseVariance	0.72	0.001	0.005		
					ClusterShade	0.62	0.017	0.032		
					LargeDependenceHighGrayLevelEmphasis	0.68	0.005	0.021		
				ngtdm	Coarseness	0.65	0.016	0.032		
					Busyness	0.64	0.012	0.032		
			Contrast		0.63	0.015	0.032			
			LLL	firstorder	10Percentile	0.65	0.002	0.017		
				glcm	InverseVariance	0.66	0.004	0.018		
				glszm	ZoneVariance	0.64	0.008	0.024		
			CTD	BIN 64	BIN 64	FO	90Percentile	0.67	0.006	0.034
						GLCM	IMC2	0.66	0.004	0.034
						GLDM	DE	0.66	0.005	0.034
	GLSZM	ZE				0.69	0.001	0.021		
	GLSZM	SAE				0.67	0.003	0.034		
LOG	Sigma: 1.0	glszm		SizeZoneNonUniformityNormalized	0.71	0.000	0.008			
				SizeZoneNonUniformityNormalized	0.71	0.000	0.008			
	Sigma: 1.5	glszm		SizeZoneNonUniformityNormalized	0.71	0.000	0.008			
				SizeZoneNonUniformityNormalized	0.75	0.000	0.001			
	Sigma: 2.0	glszm		SizeZoneNonUniformityNormalized	0.75	0.000	0.001			
				SizeZoneNonUniformityNormalized	0.75	0.000	0.001			
	Sigma: 3.0	gldm		SmallDependenceHighGrayLevelEmphasis	0.65	0.000	0.011			
				DependenceNonUniformityNormalized	0.60	0.002	0.026			
		glszm		SmallAreaHighGrayLevelEmphasis	0.63	0.002	0.026			
				SmallAreaHighGrayLevelEmphasis	0.63	0.002	0.026			
Sigma: 3.5	gldm	SmallDependenceHighGrayLevelEmphasis		0.66	0.002	0.028				
		DependenceNonUniformityNormalized		0.60	0.003	0.028				
glszm	SmallAreaHighGrayLevelEmphasis	0.65		0.002	0.028					
	SmallAreaHighGrayLevelEmphasis	0.65		0.002	0.028					
WAV	WAV	HLH		gldm	LargeDependenceHighGrayLevelEmphasis	0.67	0.002	0.038		
		HLL		firstorder	Skewness	0.67	0.002	0.045		
		LHH		firstorder	Kurtosis	0.63	0.002	0.013		
				glcm	Idn	0.68	0.001	0.013		
				glszm	GrayLevelNonUniformity	0.58	0.002	0.014		
		LLL		firstorder	90Percentile	0.68	0.004	0.034		
				glszm	SmallAreaEmphasis	0.75	0.000	0.001		
					ZoneEntropy	0.67	0.002	0.028		
			SmallAreaHighGrayLevelEmphasis		0.68	0.005	0.034			


CT	BIN64	BIN64	GLSZM	SAE	0.68	0.001	0.044
	LOG	Sigma: 1.5	glcm	InverseVariance	0.69	0.001	0.022
		Sigma: 2.0	glcm	InverseVariance	0.71	0.000	0.003
	WAV	LHL	firstorder	Skewness	0.68	0.001	0.033

Supplemental Table 3. KRAS mutation status perdition after elimination of high correlated feature in each image sets and FDR correction were presented as AUC, p-value and q-value

Image_Set	Preporcessed	Filter	Type	Features	AUC	p.value	q-value
CTD	LOG	Sigma: 3.0	glszm	LowGrayLevelZoneEmphasis	0.63	0.001	0.029
		Sigma: 4.5	firstorder	Minimum	0.71	0.001	0.010
				Skewness	0.69	0.001	0.010
			gldm	LargeDependenceLowGrayLevelEmphasis	0.68	0.001	0.013
				LowGrayLevelEmphasis	0.63	0.004	0.031
		Sigma: 4	firstorder	Minimum	0.71	0.001	0.013
				Skewness	0.68	0.001	0.013
			gldm	LargeDependenceLowGrayLevelEmphasis	0.69	0.002	0.013
				LowGrayLevelEmphasis	0.65	0.001	0.013
			glszm	LowGrayLevelZoneEmphasis	0.63	0.002	0.014
		Sigma: 5.0	firstorder	Minimum	0.70	0.001	0.014
				Skewness	0.68	0.003	0.025
			gldm	LargeDependenceLowGrayLevelEmphasis	0.71	0.000	0.000
				LowGrayLevelEmphasis	0.67	0.000	0.002

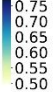

AUC, EGFR

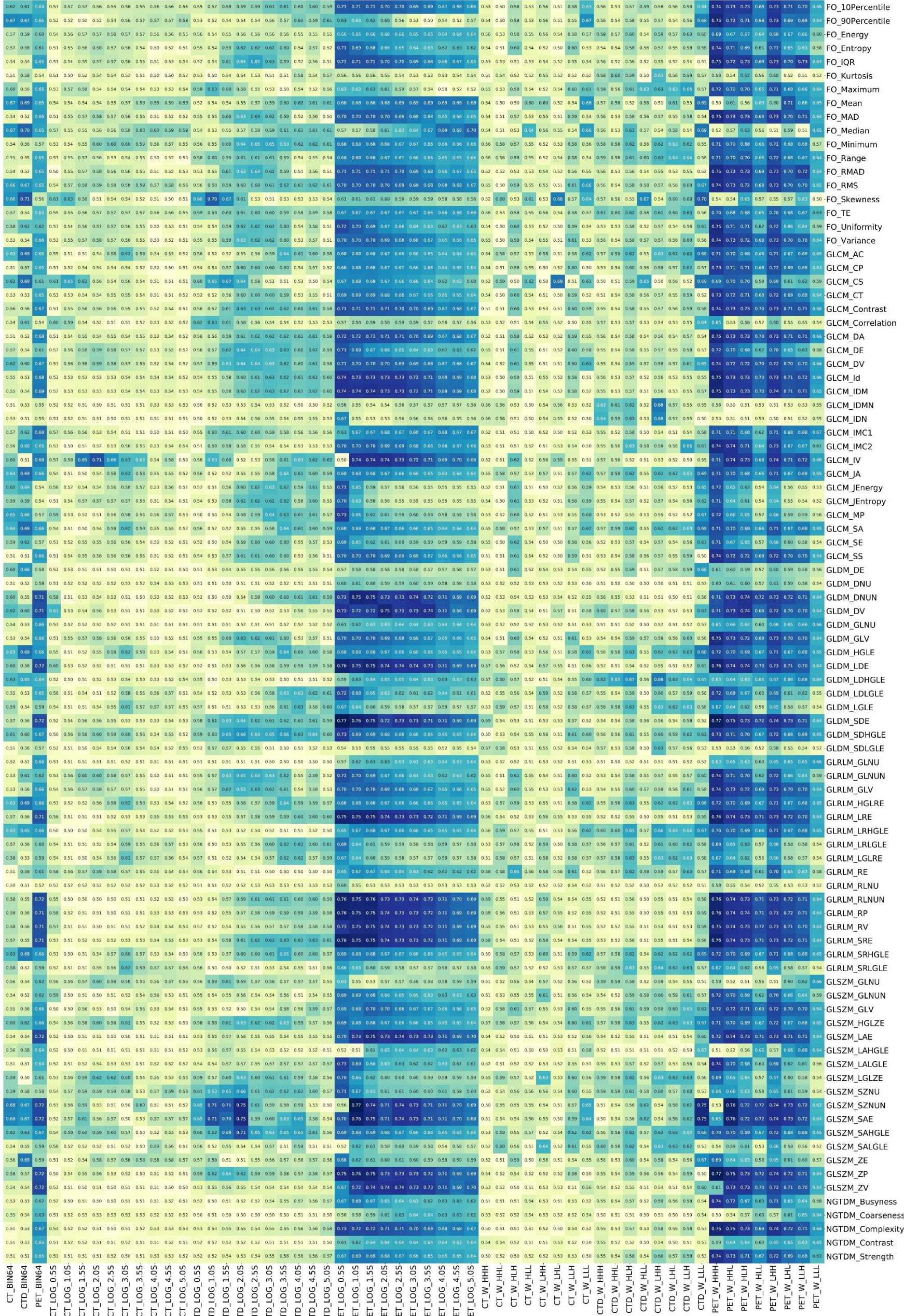


Features

Supplemental Fig 1. AUC of EGFR prediction using univariate analysis of radiomics features in different date sets




T-test, P-value, EGFR

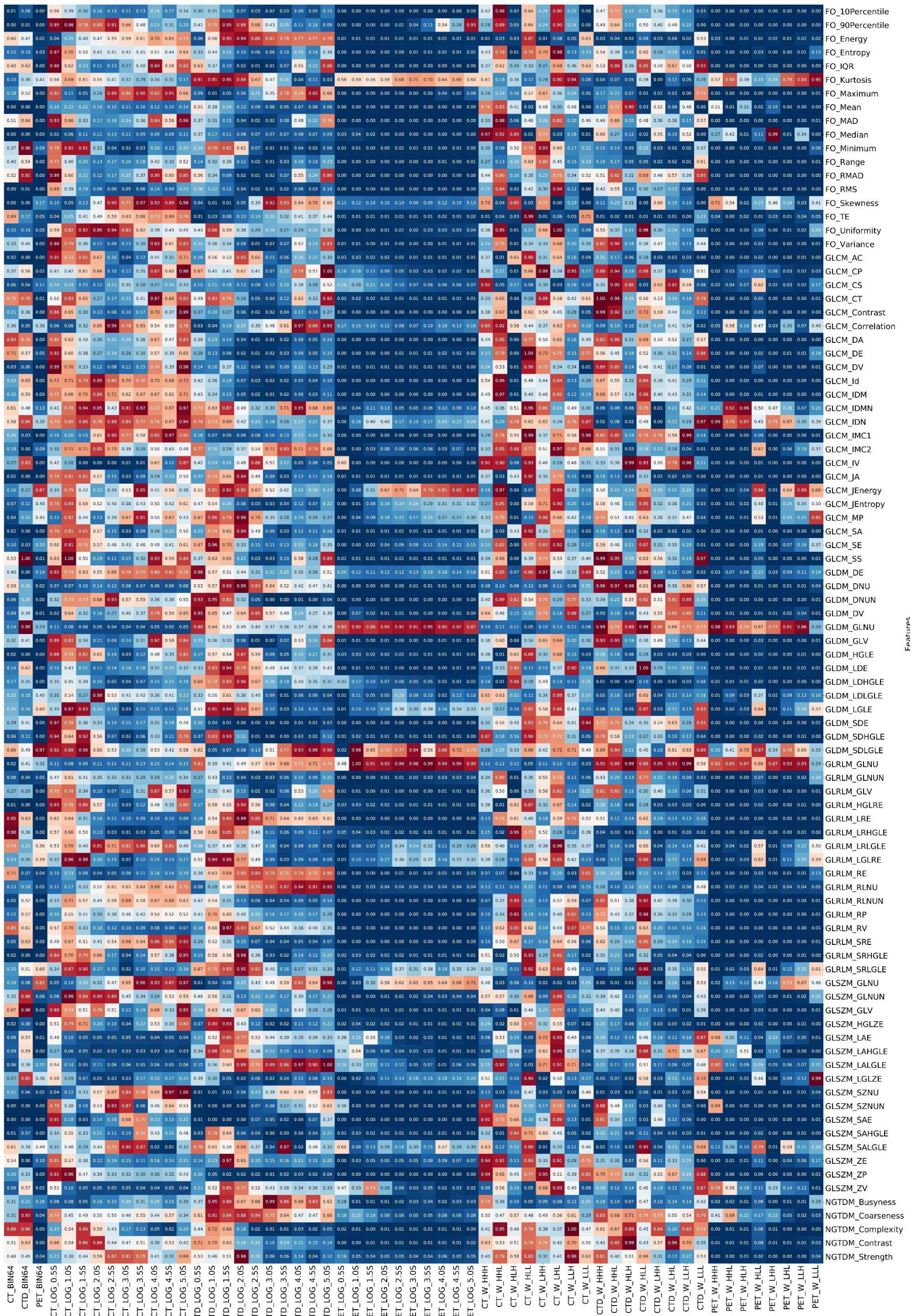


Features

Supplemental Fig 2. P-value of EGFR prediction using univariate analysis of radiomics features in different date sets

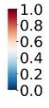


q-value, EGFR

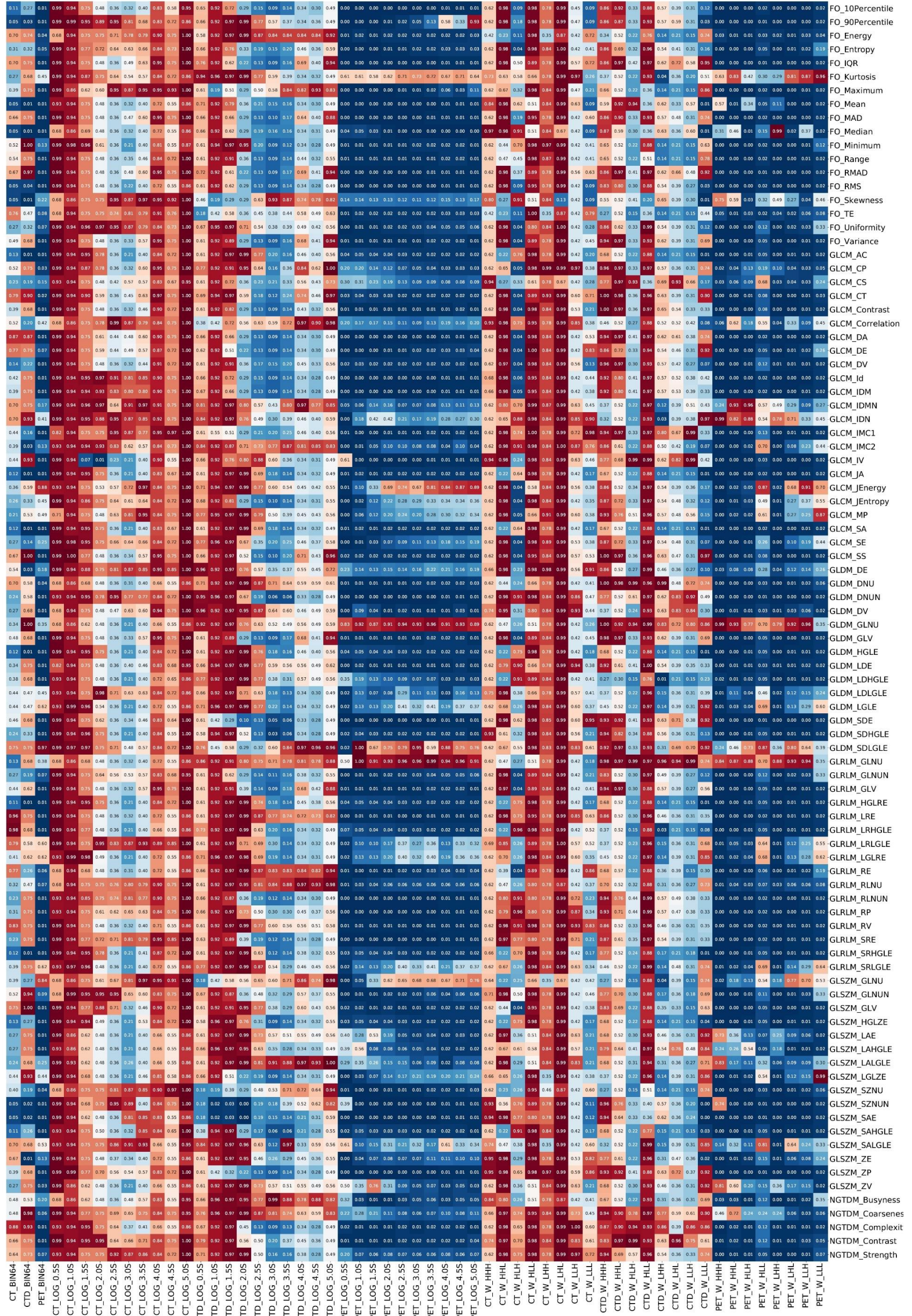


Supplemental Fig 3. Q-value of FDR correction in EGFR prediction using univariate analysis of radiomics features in different date sets




T-test, P-value, KRAS

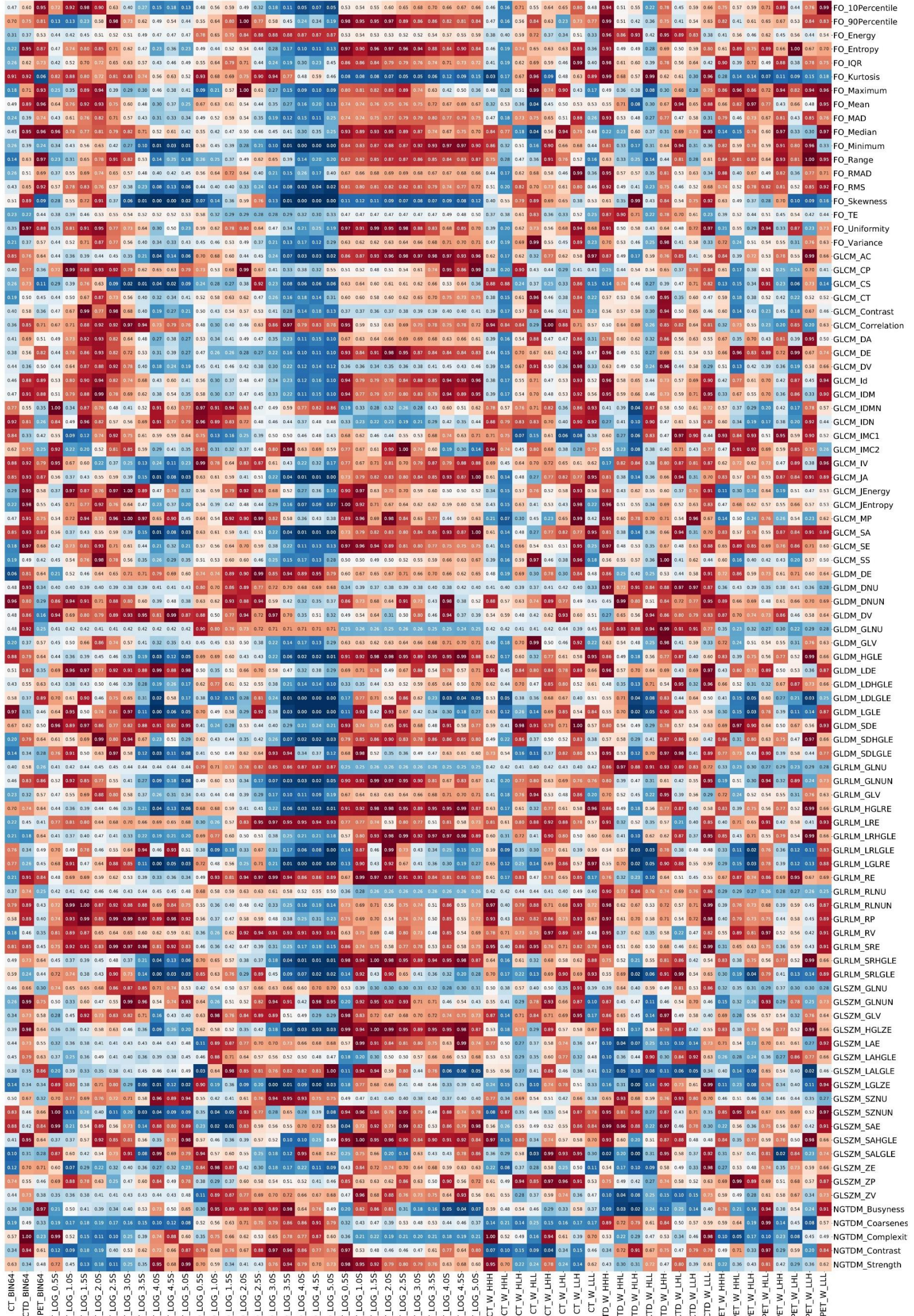

Supplemental Fig 5. P-value of KRAS prediction using univariate analysis of radiomics features in different date sets

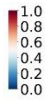


q-value, KRAS

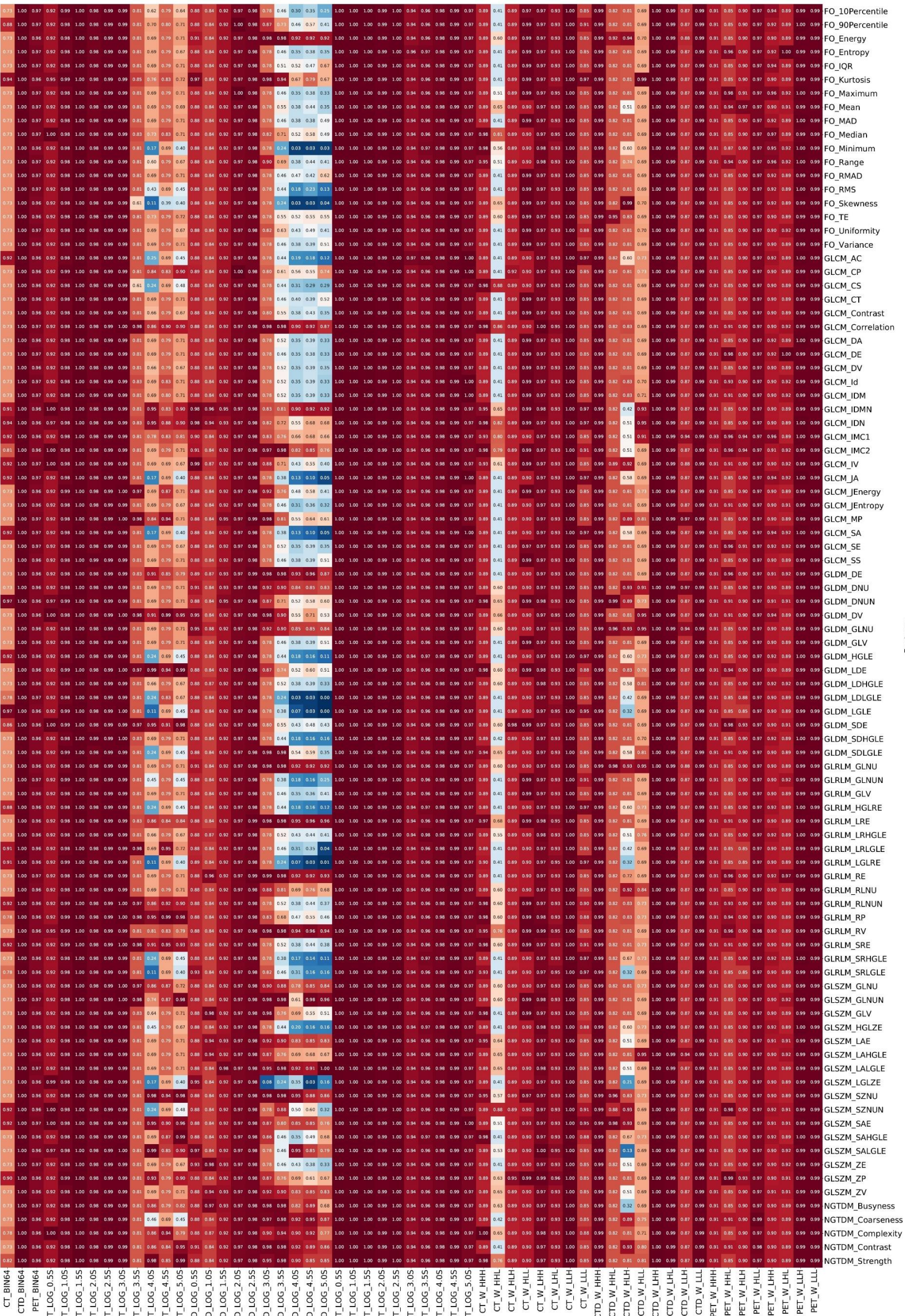

Supplemental Fig 6. Q-value of FDR correction in EGFR prediction using univariate analysis of radiomics features in different date sets

- **Shape features of EGFR and KRAS**


Details of EGFR and KRAS mutation status predation using shape features including AUC, p-value and q-value were presented in supplementary heatmap 7-9, respectively. The Surface Volume Ratio feature from PET images (AUC: 0.60, q-value: 0.28) and Flatness from CT images (AUC: 0.67, q-value:0.11) had highest perdition power in EGFR and KRAS respectively, however statistically not significant.

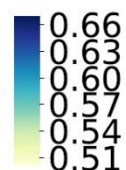

AUC, Shape Features

Feature	CT_EGFR	CTD_EGFR	PET_EGFR	CT_KRAS	CTD_KRAS	PET_KRAS
Elongation	0.53	0.54	0.52	0.62	0.59	0.60
Flatness	0.52	0.54	0.55	0.67	0.59	0.58
LeastAxis	0.52	0.53	0.58	0.51	0.51	0.54
MajorAxis	0.50	0.51	0.52	0.59	0.54	0.59
Maximum2DDiameterColumn	0.51	0.51	0.52	0.55	0.54	0.55
Maximum2DDiameterRow	0.52	0.52	0.55	0.56	0.55	0.55
Maximum2DDiameterSlice	0.52	0.53	0.54	0.54	0.54	0.55
Maximum3DDiameter	0.51	0.50	0.51	0.59	0.55	0.58
MinorAxis	0.53	0.52	0.53	0.54	0.53	0.53
Sphericity	0.55	0.57	0.57	0.58	0.54	0.60
SurfaceArea	0.52	0.52	0.55	0.53	0.53	0.56
SurfaceVolumeRatio	0.55	0.57	0.60	0.52	0.51	0.54
Volume	0.53	0.53	0.56	0.53	0.52	0.55

Supplemental Fig 7. AUC of shape features in prediction of EGFR and KRAS using univariate analysis of radiomics features in different date sets

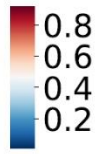


T-test, P-value, Shape Features


	CT_EGFR	CTD_EGFR	PET_EGFR	CT_KRAS	CTD_KRAS	PET_KRAS	
	0.49	0.45	0.88	0.06	0.14	0.10	Elongation
	0.56	0.35	0.28	0.01	0.10	0.19	Flatness
	0.28	0.25	0.08	0.80	0.69	0.42	LeastAxis
	0.52	0.44	0.54	0.16	0.14	0.24	MajorAxis
	0.86	0.71	0.92	0.30	0.15	0.12	Maximum2DDiameterColumn
	0.85	0.29	0.52	0.41	0.39	0.41	Maximum2DDiameterRow
	0.87	0.87	0.85	0.40	0.63	0.45	Maximum2DDiameterSlice
	0.83	0.65	0.76	0.11	0.16	0.19	Maximum3DDiameter
	0.65	0.75	0.79	0.88	0.92	0.79	MinorAxis
	0.25	0.17	0.19	0.31	0.33	0.21	Sphericity
	0.30	0.17	0.09	0.47	0.52	0.27	SurfaceArea
	0.23	0.10	0.06	0.54	0.79	0.34	SurfaceVolumeRatio
	0.12	0.05	0.05	0.43	0.47	0.30	Volume

Features

Supplemental Fig 8. P-value of shape features in prediction of EGFR and KRAS using univariate analysis of radiomics features in different date sets

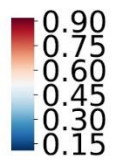

FDR, q-value, Shape Features

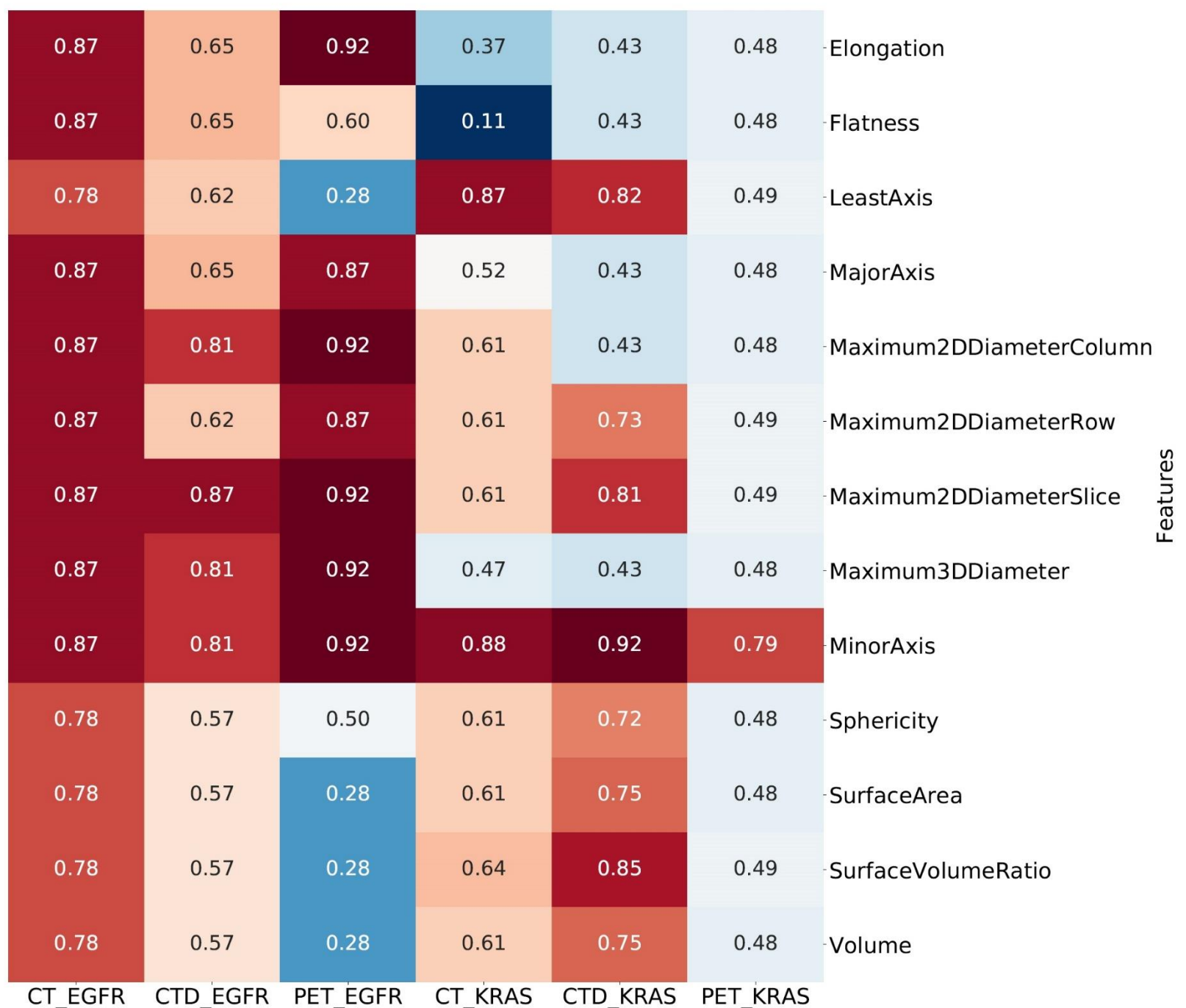

Supplemental Fig 9. Q-value of FDR correction of shape features in prediction of EGFR and KRAS using univariate analysis of radiomics features in different date sets

- Conventional PET models**

Our univariate analysis on conventional PET models for EGFR mutation status prediction showed that SUV_{peak} with AUC=0.69 (P-value = 0.0002) had the highest performance, and results for SUL_{peak} , SUV_{max} , SUL_{max} and MTV were 0.60 (P-value = 0.0003), 0.59 (P-value < 0.0001), 0.55 (P-value < 0.0001) and 0.56 (P-value = 0.0004) respectively. Our results involving conventional features for KRAS mutation status prediction were all poor: MTV (AUC: 0.55, P-value = 0.0011), SUV_{max} (AUC: 0.52, P-value < 0.0001), SUL_{max} (AUC: 0.53, P-value < 0.0001), SUV_{peak} (AUC: 0.52, P-value = 0.0013) and SUL_{peak} (AUC: 0.51, P-value = 0.0008).

Multivariate machine learning radiomic models

- EGFR**

Fig 2 shows a heatmap of EGFR mutation status prediction results in different combination of feature selection, classifier and image sets. Here, we observe a wide range of performance from 0.5 to 0.82. We further observed that the combination method of LOG preprocessed of PET image set with sigma 3.5 (PET_LOG_3.5S) with VT feature selector and SGD classifier (PET_LOG_3.5S+VT+SGD, AUC: 0.82) had the highest predictive performance, followed by PET image set which preprocessed by LHH of wavelet (PET_W_LHH) with VT_SM feature selector and MLP classifier (PET_W_LHH+ VT_SM+MLP, AUC: 0.81), LOG preprocessed of PET image set with sigma 4.0 (PET_LOG_4.0S) with SM feature selector and SGD classifier (PET_LOG_4.0S+SM+SGD, AUC: 0.81), PET_W_LHH image set with SM feature selector and SGD classifier (PET_W_LHH+SM+SGD, AUC: 0.80).

EGFR AUC in Independent Validation Set



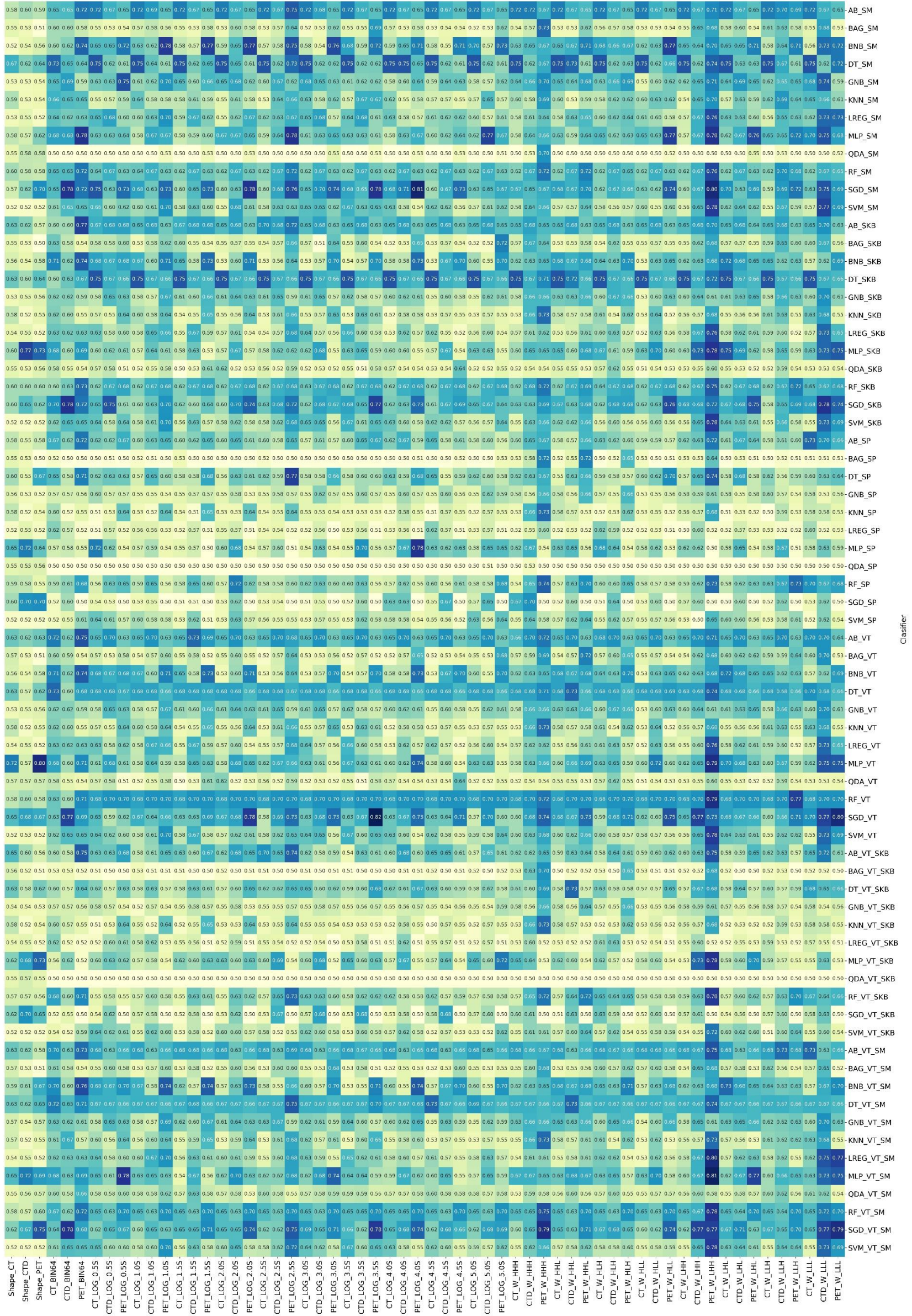



Figure 2. Heatmap depicting the predictive performance (AUC) of feature selection-classifier (rows) and image sets (columns) in prediction of EGFR mutation status of independent validation sets

Our results for EGFR mutation status prediction based on feature selection methods and image sets are depicted in supplemental Tables 4 and 5 (mean \pm SD & Min-Max). It can be deduced from that tables selection performance had a range from 0.51 to 0.82, and the combination of VT_SM feature selection with PET_W_LHH image set had the highest performance (AUC: 0.73 ± 0.074 & 0.58 - 0.81), followed by SM feature selection with PET_W_LHH image set (AUC: 0.72 ± 0.079 & 0.5 - 0.8), VT feature selection with PET_W_LHH image set (AUC: 0.71 ± 0.065 & 0.6 - 0.79) and SKB feature selection with PET_W_LHH image set (AUC: 0.7 ± 0.06 & 0.6 - 0.78). Supplemental Fig 10 delicate the box plot of EGFR mutation status prediction based on feature selection methods.

Supplemental Table 6 and 7 (mean \pm SD & Min-Max) shows our results regarding EGFR mutation status prediction based on classifier and image sets. According to these results, classifier performance has a range from 0.50 to 0.82 and combination of RF classifier with PET_W_LHH image set had the highest performance (AUC: 0.76 ± 0.024 & 0.73 - 0.79), followed by BNB classifier with PET image set which discretized into 64 Bin (PET_BIN64) (AUC: 0.75 ± 0.0096 & 0.74 - 0.76) and SVM classifier with PET_W_LHH image set (AUC: 0.75 ± 0.051 & 0.65 - 0.78) and AB classifier with PET_BIN64 image set (AUC: 0.74 ± 0.018 & 0.72 - 0.77). Supplemental figure 11 shows the box plot of EGFR mutation status prediction based on classifier method.

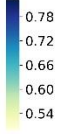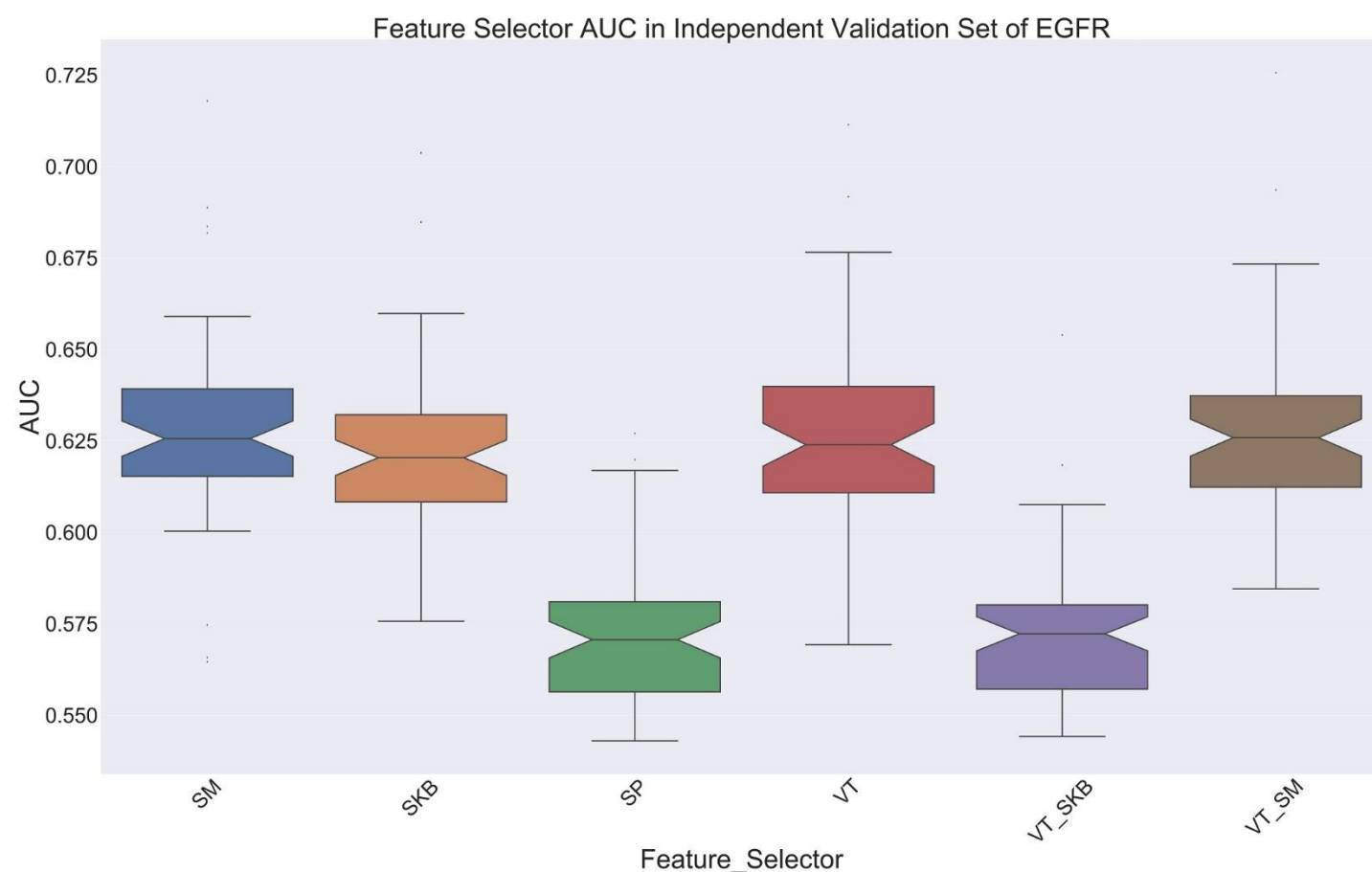

Supplemental Fig 10. AUC box plot of different feature selector in independent validation set of EGFR

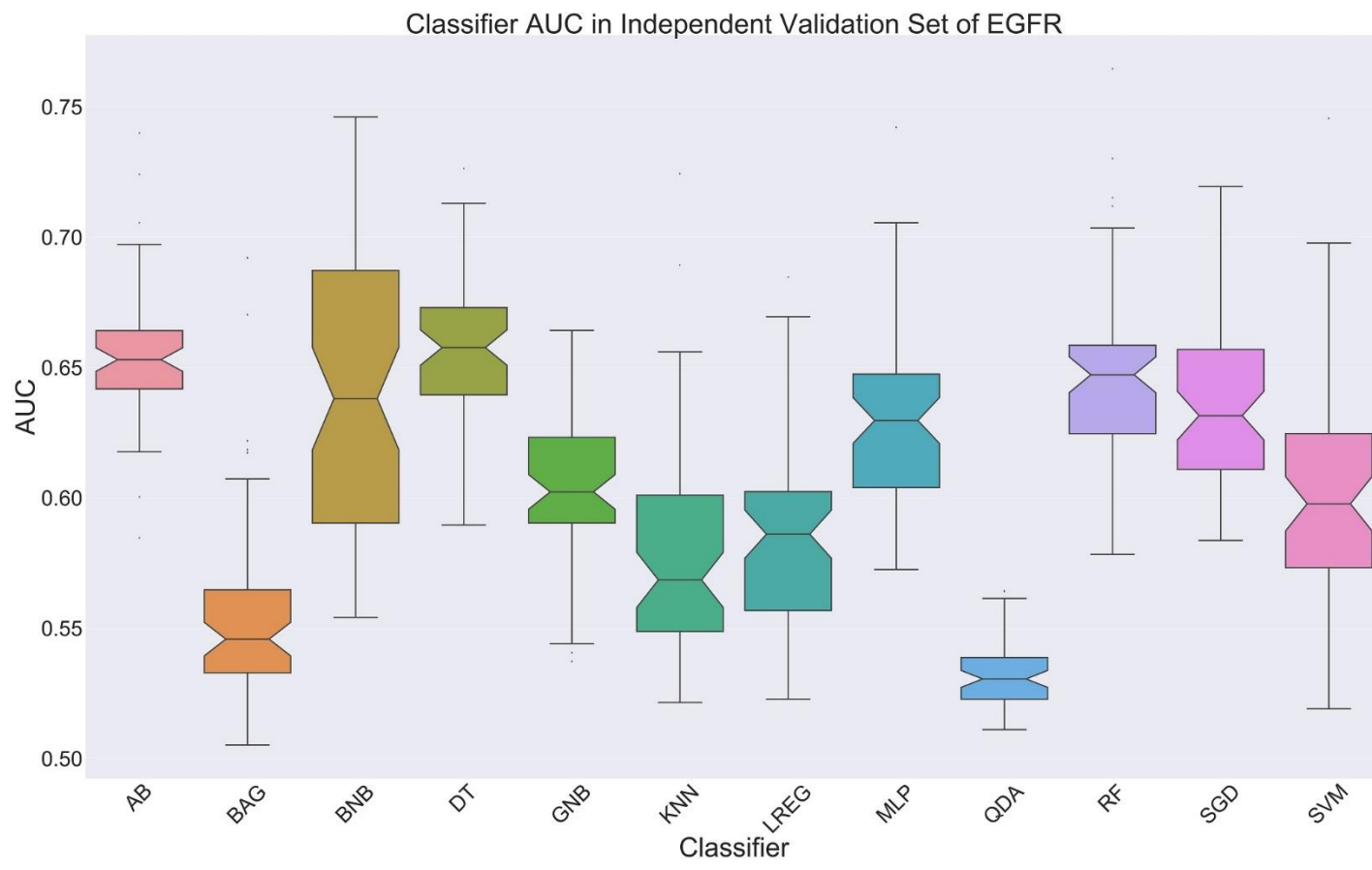

Supplemental Fig 11. AUC box plot of different classifier in independent validation set of EGFR

Supplemental Table 4. EGFR Feature Selector Min-Max in different modality, preprocessing and setting

Modality	Preprocessed	Setting	SM	SKB	SP	VT	VT_SKB	VT_SM
CT	Bin 64		0.5 - 0.73	0.58 - 0.71	0.5 - 0.67	0.58 - 0.73	0.5 - 0.68	0.6 - 0.72
	LOG	Sigma: 0.5	0.5 - 0.75	0.57 - 0.75	0.5 - 0.72	0.57 - 0.68	0.5 - 0.64	0.55 - 0.68
		Sigma: 1	0.5 - 0.75	0.52 - 0.75	0.5 - 0.6	0.52 - 0.68	0.5 - 0.58	0.53 - 0.68
		Sigma: 1.5	0.5 - 0.75	0.5 - 0.75	0.5 - 0.6	0.5 - 0.68	0.5 - 0.63	0.53 - 0.68
		sigma: 2	0.5 - 0.75	0.53 - 0.75	0.5 - 0.6	0.53 - 0.68	0.5 - 0.63	0.53 - 0.68
		sigma: 2.5	0.5 - 0.75	0.53 - 0.75	0.5 - 0.62	0.52 - 0.68	0.5 - 0.7	0.52 - 0.68
		sigma: 3	0.5 - 0.75	0.52 - 0.75	0.5 - 0.62	0.52 - 0.68	0.5 - 0.65	0.55 - 0.69
		sigma: 3.5	0.5 - 0.75	0.53 - 0.75	0.5 - 0.6	0.52 - 0.68	0.5 - 0.65	0.52 - 0.68
		sigma: 4	0.5 - 0.75	0.52 - 0.75	0.5 - 0.63	0.52 - 0.68	0.5 - 0.62	0.52 - 0.68
		sigma: 4.5	0.5 - 0.75	0.53 - 0.75	0.5 - 0.65	0.52 - 0.68	0.5 - 0.65	0.52 - 0.73
		Sigma: 5	0.5 - 0.75	0.52 - 0.75	0.5 - 0.63	0.52 - 0.68	0.5 - 0.65	0.53 - 0.69
	WAV	HHH	0.5 - 0.75	0.52 - 0.75	0.5 - 0.67	0.52 - 0.68	0.5 - 0.65	0.53 - 0.68
		HHL	0.5 - 0.75	0.53 - 0.75	0.5 - 0.63	0.54 - 0.68	0.5 - 0.62	0.53 - 0.68
		HLH	0.5 - 0.75	0.54 - 0.75	0.5 - 0.68	0.57 - 0.68	0.5 - 0.65	0.57 - 0.68
		HLL	0.5 - 0.75	0.51 - 0.75	0.5 - 0.6	0.51 - 0.68	0.5 - 0.59	0.53 - 0.68
		LHH	0.5 - 0.75	0.53 - 0.75	0.5 - 0.62	0.53 - 0.68	0.5 - 0.6	0.54 - 0.68
		LHL	0.5 - 0.75	0.55 - 0.75	0.5 - 0.61	0.55 - 0.72	0.5 - 0.6	0.55 - 0.73
		LLH	0.5 - 0.75	0.52 - 0.75	0.5 - 0.63	0.52 - 0.68	0.5 - 0.62	0.55 - 0.68
		LLL	0.5 - 0.75	0.53 - 0.75	0.5 - 0.73	0.53 - 0.7	0.5 - 0.68	0.55 - 0.73
	Shape		0.52 - 0.67	0.52 - 0.63	0.52 - 0.65	0.52 - 0.72	0.52 - 0.65	0.52 - 0.65
	CTD	Bin 64		0.5 - 0.78	0.55 - 0.78	0.5 - 0.62	0.55 - 0.77	0.5 - 0.62
LOG		Sigma: 0.5	0.5 - 0.68	0.55 - 0.75	0.5 - 0.63	0.53 - 0.7	0.5 - 0.63	0.56 - 0.67
		Sigma: 1	0.5 - 0.69	0.55 - 0.68	0.5 - 0.65	0.55 - 0.7	0.5 - 0.63	0.57 - 0.67
		Sigma: 1.5	0.5 - 0.67	0.53 - 0.68	0.5 - 0.65	0.53 - 0.73	0.5 - 0.63	0.57 - 0.68
		sigma: 2	0.5 - 0.68	0.52 - 0.7	0.5 - 0.72	0.52 - 0.7	0.5 - 0.68	0.58 - 0.7
		sigma: 2.5	0.5 - 0.68	0.52 - 0.68	0.5 - 0.6	0.52 - 0.7	0.5 - 0.69	0.55 - 0.67
		sigma: 3	0.5 - 0.7	0.51 - 0.68	0.5 - 0.63	0.53 - 0.7	0.5 - 0.68	0.53 - 0.68
		sigma: 3.5	0.5 - 0.68	0.51 - 0.68	0.5 - 0.7	0.51 - 0.7	0.5 - 0.68	0.55 - 0.67
		sigma: 4	0.5 - 0.75	0.53 - 0.68	0.5 - 0.67	0.54 - 0.7	0.5 - 0.68	0.55 - 0.68
		sigma: 4.5	0.5 - 0.67	0.54 - 0.68	0.5 - 0.67	0.53 - 0.7	0.5 - 0.68	0.53 - 0.68
		Sigma: 5	0.5 - 0.77	0.52 - 0.68	0.5 - 0.65	0.52 - 0.7	0.5 - 0.65	0.52 - 0.68
WAV		HHH	0.52 - 0.72	0.54 - 0.68	0.52 - 0.7	0.54 - 0.7	0.5 - 0.66	0.59 - 0.67
		HHL	0.5 - 0.73	0.55 - 0.72	0.5 - 0.65	0.55 - 0.73	0.5 - 0.73	0.55 - 0.73
		HLH	0.5 - 0.67	0.54 - 0.68	0.5 - 0.64	0.54 - 0.7	0.5 - 0.64	0.54 - 0.68
		HLL	0.5 - 0.67	0.54 - 0.7	0.5 - 0.62	0.54 - 0.72	0.5 - 0.6	0.57 - 0.7
		LHH	0.5 - 0.67	0.55 - 0.73	0.5 - 0.65	0.55 - 0.77	0.5 - 0.73	0.57 - 0.77
		LHL	0.5 - 0.69	0.53 - 0.69	0.5 - 0.68	0.53 - 0.7	0.5 - 0.64	0.56 - 0.71
		LLH	0.52 - 0.7	0.59 - 0.68	0.5 - 0.67	0.59 - 0.7	0.5 - 0.63	0.56 - 0.73
		LLL	0.5 - 0.77	0.53 - 0.78	0.5 - 0.7	0.53 - 0.77	0.5 - 0.72	0.62 - 0.77
Shape			0.52 - 0.62	0.52 - 0.77	0.52 - 0.72	0.52 - 0.68	0.52 - 0.7	0.52 - 0.72
PET		Bin 64		0.5 - 0.78	0.54 - 0.77	0.5 - 0.72	0.54 - 0.75	0.5 - 0.75
	LOG	Sigma: 0.5	0.5 - 0.75	0.51 - 0.68	0.5 - 0.67	0.51 - 0.7	0.5 - 0.68	0.58 - 0.78
		Sigma: 1	0.52 - 0.78	0.58 - 0.71	0.5 - 0.65	0.58 - 0.71	0.5 - 0.65	0.6 - 0.74
		Sigma: 1.5	0.53 - 0.77	0.53 - 0.73	0.5 - 0.68	0.52 - 0.73	0.5 - 0.67	0.52 - 0.74
		sigma: 2	0.53 - 0.78	0.53 - 0.74	0.5 - 0.64	0.53 - 0.78	0.5 - 0.65	0.53 - 0.74
		sigma: 2.5	0.5 - 0.78	0.59 - 0.72	0.5 - 0.77	0.59 - 0.73	0.5 - 0.74	0.55 - 0.75
		sigma: 3	0.55 - 0.76	0.52 - 0.7	0.5 - 0.66	0.52 - 0.73	0.5 - 0.62	0.55 - 0.74
		sigma: 3.5	0.52 - 0.78	0.54 - 0.77	0.5 - 0.67	0.52 - 0.82	0.5 - 0.68	0.51 - 0.78
		sigma: 4	0.52 - 0.81	0.54 - 0.73	0.5 - 0.78	0.54 - 0.74	0.5 - 0.67	0.53 - 0.74
		sigma: 4.5	0.53 - 0.73	0.52 - 0.7	0.5 - 0.63	0.52 - 0.71	0.5 - 0.61	0.55 - 0.7
		Sigma: 5	0.51 - 0.73	0.54 - 0.72	0.5 - 0.68	0.54 - 0.7	0.5 - 0.72	0.55 - 0.7
	WAV	HHH	0.64 - 0.73	0.54 - 0.73	0.5 - 0.74	0.54 - 0.74	0.5 - 0.73	0.58 - 0.79
		HHL	0.5 - 0.72	0.53 - 0.69	0.5 - 0.72	0.53 - 0.73	0.5 - 0.72	0.56 - 0.71
		HLH	0.5 - 0.69	0.55 - 0.7	0.5 - 0.66	0.55 - 0.71	0.5 - 0.66	0.54 - 0.71
		HLL	0.52 - 0.77	0.52 - 0.76	0.5 - 0.7	0.52 - 0.75	0.5 - 0.65	0.53 - 0.74
		LHH	0.5 - 0.8	0.6 - 0.78	0.5 - 0.74	0.6 - 0.79	0.5 - 0.78	0.58 - 0.81
		LHL	0.54 - 0.76	0.52 - 0.75	0.5 - 0.64	0.52 - 0.7	0.5 - 0.7	0.53 - 0.77
		LLH	0.5 - 0.72	0.52 - 0.72	0.5 - 0.73	0.52 - 0.77	0.5 - 0.7	0.54 - 0.7
		LLL	0.52 - 0.73	0.54 - 0.75	0.5 - 0.68	0.52 - 0.8	0.5 - 0.66	0.52 - 0.79
	Shape		0.51 - 0.7	0.5 - 0.73	0.5 - 0.7	0.51 - 0.8	0.51 - 0.73	0.51 - 0.75

Supplemental Table 5. EGFR Feature Selector Mean+Sd different modality, preprocessing and setting

Modality	Preprocessed	Setting	SM	SKB	SP	VT	VT_SKB	VT_SM
CT	Bin 64		0.64 ± 0.056	0.63 ± 0.042	0.58 ± 0.056	0.65 ± 0.05	0.57 ± 0.053	0.65 ± 0.042
	LOG	Sigma: 0.5	0.64 ± 0.077	0.63 ± 0.054	0.57 ± 0.07	0.63 ± 0.046	0.56 ± 0.052	0.63 ± 0.042
		Sigma: 1	0.63 ± 0.067	0.6 ± 0.062	0.56 ± 0.033	0.61 ± 0.057	0.56 ± 0.028	0.61 ± 0.047
		Sigma: 1.5	0.61 ± 0.069	0.6 ± 0.066	0.55 ± 0.032	0.6 ± 0.06	0.55 ± 0.044	0.6 ± 0.054
		sigma: 2	0.61 ± 0.077	0.6 ± 0.058	0.55 ± 0.033	0.61 ± 0.055	0.56 ± 0.043	0.61 ± 0.048
		sigma: 2.5	0.6 ± 0.072	0.6 ± 0.069	0.55 ± 0.035	0.59 ± 0.058	0.57 ± 0.057	0.6 ± 0.051
		sigma: 3	0.63 ± 0.072	0.61 ± 0.063	0.54 ± 0.036	0.61 ± 0.059	0.57 ± 0.053	0.63 ± 0.041
		sigma: 3.5	0.63 ± 0.072	0.62 ± 0.067	0.55 ± 0.036	0.62 ± 0.064	0.55 ± 0.043	0.62 ± 0.058
		sigma: 4	0.62 ± 0.069	0.6 ± 0.065	0.56 ± 0.049	0.6 ± 0.059	0.55 ± 0.045	0.61 ± 0.05
		sigma: 4.5	0.61 ± 0.07	0.59 ± 0.065	0.57 ± 0.059	0.59 ± 0.061	0.56 ± 0.053	0.61 ± 0.062
	Sigma: 5	0.61 ± 0.077	0.6 ± 0.068	0.56 ± 0.044	0.59 ± 0.056	0.55 ± 0.043	0.6 ± 0.051	
	Shape		0.57 ± 0.043	0.58 ± 0.038	0.57 ± 0.039	0.59 ± 0.057	0.58 ± 0.043	0.58 ± 0.042
	WAV	HHH	0.62 ± 0.07	0.6 ± 0.061	0.57 ± 0.051	0.6 ± 0.052	0.57 ± 0.045	0.6 ± 0.056
		HHL	0.62 ± 0.073	0.62 ± 0.062	0.56 ± 0.04	0.62 ± 0.056	0.56 ± 0.04	0.62 ± 0.055
		HLH	0.62 ± 0.069	0.62 ± 0.057	0.57 ± 0.061	0.62 ± 0.043	0.58 ± 0.051	0.63 ± 0.036
		HLL	0.61 ± 0.071	0.58 ± 0.071	0.55 ± 0.033	0.58 ± 0.063	0.55 ± 0.03	0.59 ± 0.055
		LHH	0.6 ± 0.073	0.61 ± 0.062	0.56 ± 0.041	0.6 ± 0.055	0.56 ± 0.031	0.6 ± 0.046
		LHL	0.63 ± 0.068	0.64 ± 0.072	0.55 ± 0.046	0.64 ± 0.058	0.55 ± 0.035	0.63 ± 0.057
		LLH	0.61 ± 0.069	0.61 ± 0.058	0.55 ± 0.05	0.61 ± 0.05	0.56 ± 0.04	0.62 ± 0.039
		LLL	0.63 ± 0.076	0.62 ± 0.062	0.59 ± 0.075	0.62 ± 0.062	0.58 ± 0.059	0.62 ± 0.051
	CTD	Bin 64		0.65 ± 0.065	0.62 ± 0.057	0.56 ± 0.042	0.62 ± 0.052	0.56 ± 0.035
LOG		Sigma: 0.5	0.62 ± 0.056	0.63 ± 0.058	0.57 ± 0.05	0.63 ± 0.06	0.58 ± 0.049	0.62 ± 0.035
		Sigma: 1	0.62 ± 0.051	0.62 ± 0.043	0.57 ± 0.055	0.63 ± 0.054	0.58 ± 0.047	0.62 ± 0.037
		Sigma: 1.5	0.61 ± 0.054	0.61 ± 0.053	0.56 ± 0.052	0.63 ± 0.063	0.58 ± 0.045	0.63 ± 0.038
		sigma: 2	0.64 ± 0.051	0.62 ± 0.056	0.6 ± 0.072	0.63 ± 0.06	0.59 ± 0.056	0.63 ± 0.033
		sigma: 2.5	0.63 ± 0.052	0.62 ± 0.053	0.56 ± 0.035	0.63 ± 0.061	0.58 ± 0.07	0.62 ± 0.032
		sigma: 3	0.62 ± 0.066	0.62 ± 0.064	0.57 ± 0.047	0.62 ± 0.064	0.59 ± 0.063	0.61 ± 0.045
		sigma: 3.5	0.62 ± 0.054	0.62 ± 0.051	0.59 ± 0.062	0.62 ± 0.06	0.58 ± 0.056	0.61 ± 0.035
		sigma: 4	0.62 ± 0.074	0.61 ± 0.048	0.59 ± 0.056	0.62 ± 0.054	0.59 ± 0.057	0.61 ± 0.047
		sigma: 4.5	0.61 ± 0.055	0.62 ± 0.055	0.57 ± 0.048	0.61 ± 0.061	0.58 ± 0.058	0.62 ± 0.047
Sigma: 5		0.61 ± 0.073	0.61 ± 0.059	0.58 ± 0.05	0.62 ± 0.067	0.57 ± 0.049	0.61 ± 0.051	
Shape			0.56 ± 0.035	0.58 ± 0.072	0.57 ± 0.07	0.57 ± 0.045	0.58 ± 0.063	0.58 ± 0.062
WAV		HHH	0.63 ± 0.05	0.64 ± 0.038	0.62 ± 0.056	0.64 ± 0.047	0.61 ± 0.045	0.64 ± 0.027
		HHL	0.62 ± 0.07	0.62 ± 0.056	0.58 ± 0.047	0.64 ± 0.062	0.59 ± 0.068	0.62 ± 0.053
		HLH	0.62 ± 0.053	0.63 ± 0.043	0.58 ± 0.049	0.64 ± 0.051	0.58 ± 0.052	0.63 ± 0.043
		HLL	0.6 ± 0.053	0.63 ± 0.048	0.56 ± 0.042	0.63 ± 0.059	0.56 ± 0.036	0.63 ± 0.04
		LHH	0.64 ± 0.047	0.65 ± 0.051	0.58 ± 0.052	0.65 ± 0.061	0.6 ± 0.067	0.65 ± 0.047
		LHL	0.63 ± 0.049	0.64 ± 0.055	0.59 ± 0.06	0.64 ± 0.054	0.57 ± 0.047	0.64 ± 0.043
		LLH	0.66 ± 0.048	0.64 ± 0.028	0.58 ± 0.058	0.64 ± 0.041	0.57 ± 0.041	0.64 ± 0.039
		LLL	0.69 ± 0.075	0.68 ± 0.065	0.59 ± 0.067	0.69 ± 0.063	0.6 ± 0.066	0.69 ± 0.048
PET		Bin 64		0.66 ± 0.081	0.65 ± 0.083	0.58 ± 0.087	0.65 ± 0.079	0.58 ± 0.086
	LOG	Sigma: 0.5	0.63 ± 0.075	0.63 ± 0.048	0.57 ± 0.064	0.63 ± 0.049	0.58 ± 0.064	0.65 ± 0.057
		Sigma: 1	0.65 ± 0.074	0.65 ± 0.041	0.56 ± 0.054	0.65 ± 0.044	0.56 ± 0.052	0.66 ± 0.043
		Sigma: 1.5	0.63 ± 0.065	0.63 ± 0.062	0.56 ± 0.071	0.64 ± 0.06	0.57 ± 0.067	0.64 ± 0.067
		sigma: 2	0.64 ± 0.074	0.62 ± 0.069	0.56 ± 0.055	0.64 ± 0.074	0.57 ± 0.066	0.64 ± 0.064
		sigma: 2.5	0.68 ± 0.08	0.66 ± 0.044	0.57 ± 0.086	0.66 ± 0.041	0.59 ± 0.091	0.67 ± 0.062
		sigma: 3	0.64 ± 0.059	0.61 ± 0.059	0.55 ± 0.05	0.62 ± 0.069	0.55 ± 0.045	0.65 ± 0.061
		sigma: 3.5	0.65 ± 0.065	0.63 ± 0.065	0.55 ± 0.053	0.63 ± 0.079	0.57 ± 0.065	0.64 ± 0.075
		sigma: 4	0.63 ± 0.078	0.62 ± 0.066	0.57 ± 0.088	0.64 ± 0.073	0.56 ± 0.051	0.64 ± 0.071
		sigma: 4.5	0.62 ± 0.065	0.61 ± 0.065	0.55 ± 0.049	0.61 ± 0.072	0.54 ± 0.04	0.61 ± 0.059
	Sigma: 5	0.62 ± 0.054	0.62 ± 0.062	0.57 ± 0.063	0.62 ± 0.06	0.57 ± 0.067	0.62 ± 0.057	
	Shape		0.57 ± 0.057	0.58 ± 0.064	0.57 ± 0.069	0.59 ± 0.081	0.57 ± 0.069	0.6 ± 0.076
	WAV	HHH	0.68 ± 0.028	0.66 ± 0.051	0.63 ± 0.088	0.68 ± 0.056	0.62 ± 0.092	0.67 ± 0.059
		HHL	0.64 ± 0.061	0.64 ± 0.052	0.6 ± 0.083	0.65 ± 0.059	0.57 ± 0.074	0.64 ± 0.049
		HLH	0.62 ± 0.051	0.62 ± 0.054	0.58 ± 0.061	0.63 ± 0.054	0.58 ± 0.062	0.63 ± 0.055
		HLL	0.64 ± 0.083	0.61 ± 0.071	0.55 ± 0.06	0.61 ± 0.076	0.55 ± 0.052	0.63 ± 0.065
		LHH	0.72 ± 0.079	0.7 ± 0.06	0.62 ± 0.098	0.71 ± 0.065	0.65 ± 0.11	0.73 ± 0.074
		LHL	0.64 ± 0.064	0.62 ± 0.063	0.55 ± 0.048	0.62 ± 0.049	0.57 ± 0.068	0.63 ± 0.067
		LLH	0.64 ± 0.066	0.61 ± 0.06	0.56 ± 0.07	0.62 ± 0.071	0.56 ± 0.064	0.63 ± 0.055
		LLL	0.65 ± 0.072	0.64 ± 0.07	0.57 ± 0.066	0.65 ± 0.084	0.56 ± 0.061	0.66 ± 0.091

Supplemental Table 6. EGFR Classifier Min-Max different modality, preprocessing and setting

Modality	Preprocessed	Setting	AB	BAG	BNB	DT	GNB	KNN	LREG	MLP	QDA	RF	SGD	SVM
CT	Bin 64		0.6 - 0.72	0.52 - 0.63	0.6 - 0.71	0.6 - 0.73	0.57 - 0.65	0.6 - 0.66	0.62 - 0.64	0.56 - 0.68	0.5 - 0.6	0.59 - 0.68	0.52 - 0.7	0.52 - 0.62
	LOG	Sigma: 0.5	0.62 - 0.72	0.5 - 0.58	0.65 - 0.68	0.62 - 0.75	0.58 - 0.63	0.51 - 0.6	0.51 - 0.65	0.6 - 0.72	0.5 - 0.58	0.55 - 0.68	0.54 - 0.75	0.64 - 0.66
		Sigma: 1	0.58 - 0.72	0.52 - 0.57	0.63 - 0.67	0.57 - 0.75	0.52 - 0.61	0.53 - 0.64	0.56 - 0.6	0.57 - 0.63	0.5 - 0.58	0.57 - 0.68	0.53 - 0.68	0.55 - 0.6
		Sigma: 1.5	0.6 - 0.72	0.51 - 0.6	0.58 - 0.65	0.58 - 0.75	0.55 - 0.65	0.52 - 0.58	0.52 - 0.59	0.54 - 0.62	0.5 - 0.57	0.55 - 0.68	0.51 - 0.65	0.53 - 0.58
		sigma: 2	0.6 - 0.72	0.5 - 0.55	0.53 - 0.59	0.56 - 0.75	0.55 - 0.65	0.53 - 0.55	0.52 - 0.6	0.57 - 0.67	0.5 - 0.62	0.55 - 0.68	0.53 - 0.67	0.53 - 0.6
		sigma: 2.5	0.6 - 0.72	0.51 - 0.54	0.53 - 0.58	0.62 - 0.75	0.55 - 0.63	0.53 - 0.54	0.54 - 0.62	0.57 - 0.62	0.5 - 0.6	0.57 - 0.68	0.53 - 0.63	0.56 - 0.58
		sigma: 3	0.57 - 0.72	0.5 - 0.57	0.53 - 0.6	0.58 - 0.75	0.55 - 0.68	0.55 - 0.63	0.52 - 0.65	0.54 - 0.63	0.5 - 0.57	0.62 - 0.68	0.51 - 0.69	0.53 - 0.65
		sigma: 3.5	0.54 - 0.72	0.5 - 0.55	0.53 - 0.68	0.58 - 0.75	0.55 - 0.64	0.53 - 0.57	0.53 - 0.66	0.55 - 0.67	0.5 - 0.59	0.58 - 0.68	0.52 - 0.68	0.53 - 0.67
		sigma: 4	0.6 - 0.72	0.5 - 0.57	0.58 - 0.6	0.62 - 0.75	0.55 - 0.6	0.52 - 0.62	0.51 - 0.63	0.57 - 0.6	0.5 - 0.57	0.57 - 0.68	0.53 - 0.68	0.57 - 0.63
		sigma: 4.5	0.65 - 0.72	0.5 - 0.53	0.53 - 0.58	0.63 - 0.75	0.55 - 0.59	0.5 - 0.57	0.55 - 0.62	0.57 - 0.63	0.5 - 0.6	0.6 - 0.68	0.55 - 0.63	0.52 - 0.62
	Sigma: 5	0.57 - 0.72	0.5 - 0.55	0.6 - 0.7	0.58 - 0.75	0.55 - 0.58	0.52 - 0.57	0.52 - 0.6	0.57 - 0.65	0.5 - 0.58	0.57 - 0.68	0.57 - 0.65	0.53 - 0.63	
	Shape		0.58 - 0.65	0.55 - 0.57	0.52 - 0.59	0.6 - 0.67	0.53 - 0.57	0.57 - 0.59	0.52 - 0.54	0.58 - 0.72	0.55 - 0.57	0.57 - 0.6	0.57 - 0.65	0.52 - 0.52
	WAV	HHH	0.6 - 0.72	0.53 - 0.57	0.62 - 0.63	0.61 - 0.75	0.58 - 0.64	0.53 - 0.6	0.52 - 0.58	0.58 - 0.67	0.5 - 0.53	0.54 - 0.68	0.6 - 0.67	0.55 - 0.6
		HHL	0.58 - 0.72	0.5 - 0.54	0.65 - 0.68	0.55 - 0.75	0.58 - 0.65	0.58 - 0.6	0.52 - 0.58	0.6 - 0.65	0.5 - 0.56	0.57 - 0.68	0.51 - 0.68	0.57 - 0.62
		HLH	0.58 - 0.72	0.5 - 0.6	0.64 - 0.68	0.59 - 0.75	0.57 - 0.63	0.52 - 0.61	0.6 - 0.62	0.62 - 0.68	0.5 - 0.57	0.6 - 0.68	0.51 - 0.67	0.55 - 0.6
		HLL	0.59 - 0.72	0.52 - 0.57	0.53 - 0.62	0.57 - 0.75	0.53 - 0.55	0.53 - 0.6	0.52 - 0.61	0.58 - 0.65	0.5 - 0.55	0.58 - 0.68	0.52 - 0.63	0.55 - 0.6
		LHH	0.6 - 0.72	0.52 - 0.55	0.62 - 0.65	0.57 - 0.75	0.58 - 0.62	0.54 - 0.56	0.5 - 0.58	0.53 - 0.62	0.5 - 0.56	0.59 - 0.68	0.57 - 0.68	0.53 - 0.56
		LHL	0.58 - 0.72	0.5 - 0.6	0.65 - 0.73	0.58 - 0.75	0.57 - 0.64	0.51 - 0.57	0.52 - 0.63	0.58 - 0.75	0.5 - 0.55	0.57 - 0.68	0.5 - 0.7	0.6 - 0.64
		LLH	0.58 - 0.72	0.5 - 0.62	0.58 - 0.65	0.6 - 0.75	0.58 - 0.62	0.5 - 0.62	0.53 - 0.59	0.58 - 0.65	0.5 - 0.6	0.57 - 0.68	0.52 - 0.67	0.51 - 0.55
		LLL	0.65 - 0.73	0.51 - 0.6	0.56 - 0.57	0.6 - 0.75	0.58 - 0.68	0.58 - 0.65	0.57 - 0.62	0.55 - 0.7	0.5 - 0.61	0.62 - 0.7	0.53 - 0.7	0.52 - 0.57
CTD	Bin 65		0.58 - 0.65	0.5 - 0.6	0.6 - 0.62	0.57 - 0.65	0.57 - 0.69	0.52 - 0.67	0.52 - 0.63	0.58 - 0.68	0.5 - 0.58	0.6 - 0.65	0.55 - 0.78	0.52 - 0.65
	LOG	Sigma: 0.5	0.62 - 0.7	0.5 - 0.6	0.65 - 0.67	0.57 - 0.68	0.57 - 0.65	0.53 - 0.57	0.57 - 0.68	0.61 - 0.68	0.5 - 0.58	0.58 - 0.7	0.53 - 0.75	0.58 - 0.62
		Sigma: 1	0.61 - 0.7	0.51 - 0.58	0.58 - 0.62	0.63 - 0.68	0.55 - 0.62	0.52 - 0.6	0.56 - 0.67	0.58 - 0.67	0.5 - 0.58	0.6 - 0.7	0.55 - 0.65	0.58 - 0.62
		Sigma: 1.5	0.6 - 0.73	0.5 - 0.58	0.57 - 0.58	0.58 - 0.68	0.57 - 0.6	0.51 - 0.61	0.56 - 0.67	0.57 - 0.67	0.5 - 0.63	0.63 - 0.7	0.51 - 0.65	0.58 - 0.67
		sigma: 2	0.63 - 0.7	0.5 - 0.6	0.6 - 0.65	0.6 - 0.68	0.58 - 0.68	0.53 - 0.61	0.57 - 0.62	0.6 - 0.7	0.5 - 0.58	0.63 - 0.72	0.62 - 0.7	0.6 - 0.68
		sigma: 2.5	0.58 - 0.7	0.5 - 0.6	0.55 - 0.58	0.59 - 0.68	0.55 - 0.67	0.55 - 0.64	0.54 - 0.63	0.6 - 0.69	0.5 - 0.58	0.58 - 0.7	0.54 - 0.69	0.52 - 0.63
		sigma: 3	0.58 - 0.7	0.5 - 0.53	0.54 - 0.57	0.58 - 0.68	0.58 - 0.65	0.55 - 0.58	0.54 - 0.68	0.61 - 0.68	0.5 - 0.58	0.63 - 0.7	0.55 - 0.7	0.57 - 0.65
		sigma: 3.5	0.63 - 0.7	0.5 - 0.6	0.55 - 0.59	0.6 - 0.68	0.56 - 0.62	0.53 - 0.67	0.56 - 0.68	0.6 - 0.7	0.5 - 0.59	0.62 - 0.7	0.6 - 0.68	0.6 - 0.63
		sigma: 4	0.63 - 0.7	0.5 - 0.57	0.55 - 0.65	0.59 - 0.75	0.56 - 0.62	0.55 - 0.58	0.56 - 0.62	0.6 - 0.67	0.5 - 0.57	0.58 - 0.7	0.58 - 0.71	0.58 - 0.62
		sigma: 4.5	0.58 - 0.7	0.5 - 0.57	0.55 - 0.67	0.55 - 0.68	0.54 - 0.64	0.57 - 0.57	0.55 - 0.62	0.6 - 0.67	0.5 - 0.58	0.56 - 0.7	0.64 - 0.68	0.56 - 0.62
	Sigma: 5	0.6 - 0.7	0.5 - 0.53	0.55 - 0.57	0.6 - 0.68	0.54 - 0.62	0.57 - 0.65	0.57 - 0.63	0.6 - 0.77	0.5 - 0.57	0.58 - 0.7	0.6 - 0.7	0.52 - 0.61	
	Shape		0.55 - 0.62	0.52 - 0.53	0.54 - 0.61	0.55 - 0.62	0.53 - 0.55	0.52 - 0.53	0.55 - 0.55	0.57 - 0.77	0.53 - 0.58	0.57 - 0.6	0.62 - 0.7	0.52 - 0.53
	WAV	HHH	0.62 - 0.72	0.57 - 0.67	0.63 - 0.65	0.6 - 0.68	0.56 - 0.66	0.58 - 0.66	0.55 - 0.62	0.63 - 0.67	0.5 - 0.59	0.65 - 0.7	0.61 - 0.7	0.61 - 0.64
		HHL	0.57 - 0.7	0.52 - 0.57	0.67 - 0.67	0.63 - 0.73	0.56 - 0.64	0.53 - 0.57	0.52 - 0.63	0.59 - 0.68	0.5 - 0.6	0.63 - 0.7	0.6 - 0.67	0.57 - 0.62
		HLH	0.62 - 0.7	0.5 - 0.62	0.63 - 0.66	0.58 - 0.68	0.55 - 0.67	0.53 - 0.62	0.59 - 0.63	0.57 - 0.65	0.5 - 0.62	0.6 - 0.7	0.59 - 0.68	0.56 - 0.62
		HLL	0.59 - 0.7	0.5 - 0.6	0.63 - 0.63	0.57 - 0.68	0.55 - 0.6	0.55 - 0.64	0.53 - 0.63	0.6 - 0.72	0.5 - 0.57	0.57 - 0.7	0.6 - 0.63	0.55 - 0.6
		LHH	0.63 - 0.7	0.52 - 0.65	0.63 - 0.64	0.62 - 0.68	0.59 - 0.65	0.57 - 0.65	0.6 - 0.67	0.62 - 0.73	0.5 - 0.6	0.62 - 0.7	0.6 - 0.77	0.5 - 0.65
		LHL	0.59 - 0.7	0.5 - 0.62	0.65 - 0.68	0.63 - 0.68	0.55 - 0.69	0.53 - 0.63	0.55 - 0.63	0.6 - 0.69	0.5 - 0.58	0.6 - 0.7	0.56 - 0.71	0.6 - 0.64
		LLH	0.61 - 0.73	0.52 - 0.65	0.62 - 0.64	0.56 - 0.68	0.57 - 0.66	0.59 - 0.69	0.52 - 0.65	0.57 - 0.67	0.5 - 0.62	0.63 - 0.7	0.6 - 0.69	0.58 - 0.67
		LLL	0.63 - 0.72	0.51 - 0.7	0.62 - 0.73	0.62 - 0.68	0.53 - 0.74	0.58 - 0.68	0.52 - 0.75	0.63 - 0.75	0.5 - 0.62	0.64 - 0.72	0.62 - 0.78	0.6 - 0.77
PET	Bin 66		0.72 - 0.77	0.5 - 0.58	0.74 - 0.76	0.64 - 0.71	0.56 - 0.61	0.55 - 0.65	0.52 - 0.63	0.52 - 0.78	0.5 - 0.66	0.68 - 0.73	0.5 - 0.72	0.59 - 0.65
	LOG	Sigma: 0.5	0.63 - 0.68	0.5 - 0.6	0.68 - 0.72	0.61 - 0.67	0.54 - 0.75	0.59 - 0.64	0.52 - 0.62	0.54 - 0.78	0.5 - 0.6	0.55 - 0.7	0.5 - 0.73	0.57 - 0.64
		Sigma: 1	0.63 - 0.66	0.5 - 0.62	0.71 - 0.78	0.57 - 0.66	0.54 - 0.7	0.58 - 0.64	0.53 - 0.7	0.54 - 0.67	0.5 - 0.62	0.56 - 0.7	0.5 - 0.73	0.58 - 0.7
		Sigma: 1.5	0.65 - 0.69	0.5 - 0.58	0.73 - 0.77	0.61 - 0.68	0.57 - 0.67	0.59 - 0.65	0.51 - 0.62	0.5 - 0.63	0.5 - 0.61	0.6 - 0.7	0.5 - 0.73	0.51 - 0.6
		sigma: 2	0.61 - 0.66	0.5 - 0.58	0.71 - 0.77	0.61 - 0.66	0.53 - 0.63	0.58 - 0.64	0.51 - 0.67	0.54 - 0.65	0.5 - 0.53	0.62 - 0.7	0.5 - 0.78	0.57 - 0.61
		sigma: 2.5	0.65 - 0.75	0.5 - 0.66	0.64 - 0.75	0.65 - 0.77	0.57 - 0.62	0.64 - 0.68	0.52 - 0.68	0.51 - 0.78	0.5 - 0.59	0.6 - 0.73	0.5 - 0.76	0.55 - 0.72
		sigma: 3	0.57 - 0.66	0.5 - 0.68	0.7 - 0.76	0.62 - 0.66	0.57 - 0.63	0.54 - 0.65	0.5 - 0.57	0.52 - 0.74	0.5 - 0.59	0.6 - 0.7	0.5 - 0.74	0.55 - 0.62
		sigma: 3.5	0.61 - 0.67	0.5 - 0.69	0.7 - 0.72	0.6 - 0.7	0.57 - 0.62	0.53 - 0.67	0.51 - 0.61	0.56 - 0.63	0.5 - 0.58	0.56 - 0.7	0.5 - 0.82	0.58 - 0.65
		sigma: 4	0.59 - 0.66	0.51 - 0.65	0.71 - 0.74	0.65 - 0.68	0.57 - 0.65	0.55 - 0.6	0.51 - 0.57	0.55 - 0.78	0.5 - 0.58	0.56 - 0.7	0.5 - 0.81	0.54 - 0.58
		sigma: 4.5	0.61 - 0.66	0.51 - 0.6	0.7 - 0.71	0.59 - 0.66	0.56 - 0.63	0.55 - 0.55	0.51 - 0.62	0.54 - 0.64	0.5 - 0.64	0.59 - 0.7	0.5 - 0.73	0.52 - 0.57
	Sigma: 5	0.56 - 0.66	0.52 - 0.72	0.7 - 0.73	0.58 - 0.66	0.59 - 0.62	0.55 - 0.57	0.51 - 0.61	0.55 - 0.72	0.5 - 0.58	0.58 - 0.7	0.5 - 0.69	0.62 - 0.64	
	Shape		0.56 - 0.63	0.5 - 0.51	0.56 - 0.67	0.62 - 0.67	0.52 - 0.57	0.54 - 0.55	0.52 - 0.52	0.62 - 0.8	0.55 - 0.58	0.55 - 0.6	0.62 - 0.75	0.52 - 0.52
	WAV	HHH	0.65 - 0.72	0.64 - 0.73	0.65 - 0.67	0.66 - 0.71	0.66 - 0.7	0.69 - 0.73	0.52 - 0.64	0.52 - 0.66	0.5 - 0.7	0.72 - 0.75	0.5 - 0.79	0.58 - 0.68
		HHL	0.63 - 0.66	0.52 - 0.72	0.68 - 0.71	0.57 - 0.66	0.64 - 0.68	0.53 - 0.61	0.52 - 0.65	0.54 - 0.69	0.5 - 0.56	0.69 - 0.72	0.5 - 0.73	0.62 - 0.66
		HLH	0.6 - 0.66	0.55 - 0.65	0.67 - 0.71	0.57 - 0.66	0.66 - 0.69	0.62 - 0.62	0.52 - 0.64	0.52 - 0.63	0.5 - 0.55	0.65 - 0.7	0.5 - 0.71	0.54 - 0.6
		HLL	0.57 - 0.66	0.51 - 0.58	0.65 - 0.77	0.65 - 0.7	0.56 - 0.66	0.52 - 0.62	0.51 - 0.58	0.53 - 0.77	0.5 - 0.57	0.58 - 0.7	0.5 - 0.76	0.56 - 0.6
		LHH	0.7 - 0.75	0.64 - 0.68	0.68 - 0.7	0.68 - 0.74	0.61 - 0.71	0.68 - 0.73	0.52 - 0.8	0.5 - 0.81	0.5 - 0.6	0.73 - 0.79	0.5 - 0.8	0.65 - 0.78
		LHL	0.63 - 0.66	0.5 - 0.62	0.65 - 0.71	0.57 - 0.66	0.56 - 0.66	0.52 - 0.59	0.53 - 0.63	0.54 - 0.77	0.5 - 0.57	0.62 - 0.7	0.5 - 0.75	0.56 - 0.62
		LLH	0.57 - 0.69	0.5 - 0.64	0.63 - 0.71	0.59 - 0.66	0.54 - 0.65	0.53 - 0.64	0.52 - 0.63	0.51 - 0.72	0.5 - 0.56	0.68 - 0.77	0.5 - 0.72	0.58 - 0.64
		LLL	0.61 - 0.66	0.5 - 0.56										

Supplemental Table 7. EGFR Classifier Mean+Sd different modality, preprocessing and settin

Modality	Preprocessed	Setting	AB	BAG	BNB	DT	GNB	KNN	LREG	MLP	QDA	RF	SGD	SVM	
CT	Bin 64		0.66 ± 0.049	0.58 ± 0.048	0.68 ± 0.054	0.67 ± 0.064	0.61 ± 0.035	0.62 ± 0.024	0.63 ± 0.01	0.64 ± 0.063	0.54 ± 0.049	0.64 ± 0.036	0.61 ± 0.073	0.59 ± 0.042	
	LOG	Sigma: 0.5	0.66 ± 0.034	0.54 ± 0.032	0.67 ± 0.017	0.68 ± 0.06	0.59 ± 0.018	0.55 ± 0.034	0.6 ± 0.066	0.64 ± 0.044	0.54 ± 0.04	0.62 ± 0.052	0.62 ± 0.078	0.64 ± 0.0082	
		Sigma: 1	0.65 ± 0.05	0.54 ± 0.016	0.66 ± 0.017	0.67 ± 0.079	0.57 ± 0.031	0.58 ± 0.038	0.58 ± 0.016	0.58 ± 0.026	0.52 ± 0.032	0.62 ± 0.042	0.61 ± 0.057	0.58 ± 0.018	
		Sigma: 1.5	0.66 ± 0.04	0.55 ± 0.03	0.62 ± 0.032	0.67 ± 0.075	0.6 ± 0.04	0.54 ± 0.022	0.56 ± 0.024	0.57 ± 0.027	0.51 ± 0.027	0.62 ± 0.047	0.59 ± 0.061	0.56 ± 0.024	
		sigma: 2	0.65 ± 0.043	0.53 ± 0.021	0.56 ± 0.027	0.66 ± 0.084	0.6 ± 0.039	0.54 ± 0.0091	0.56 ± 0.025	0.61 ± 0.034	0.55 ± 0.058	0.61 ± 0.05	0.6 ± 0.056	0.57 ± 0.028	
		sigma: 2.5	0.68 ± 0.043	0.52 ± 0.013	0.55 ± 0.025	0.68 ± 0.06	0.59 ± 0.035	0.53 ± 0.0017	0.56 ± 0.03	0.6 ± 0.02	0.54 ± 0.042	0.62 ± 0.043	0.58 ± 0.042	0.58 ± 0.011	
		sigma: 3	0.65 ± 0.052	0.53 ± 0.026	0.56 ± 0.034	0.68 ± 0.064	0.6 ± 0.049	0.58 ± 0.039	0.6 ± 0.063	0.6 ± 0.031	0.52 ± 0.029	0.64 ± 0.025	0.61 ± 0.062	0.62 ± 0.044	
		sigma: 3.5	0.64 ± 0.063	0.52 ± 0.021	0.57 ± 0.073	0.67 ± 0.073	0.6 ± 0.038	0.54 ± 0.014	0.61 ± 0.062	0.63 ± 0.04	0.53 ± 0.036	0.62 ± 0.036	0.62 ± 0.073	0.62 ± 0.069	
		sigma: 4	0.65 ± 0.041	0.52 ± 0.025	0.59 ± 0.0088	0.68 ± 0.056	0.58 ± 0.025	0.54 ± 0.04	0.55 ± 0.045	0.58 ± 0.015	0.53 ± 0.037	0.63 ± 0.039	0.63 ± 0.05	0.61 ± 0.031	
		sigma: 4.5	0.67 ± 0.028	0.52 ± 0.011	0.55 ± 0.025	0.7 ± 0.052	0.56 ± 0.014	0.53 ± 0.024	0.6 ± 0.026	0.59 ± 0.025	0.53 ± 0.039	0.63 ± 0.031	0.6 ± 0.029	0.59 ± 0.044	
	Sigma: 5	0.64 ± 0.057	0.53 ± 0.02	0.62 ± 0.05	0.68 ± 0.068	0.56 ± 0.014	0.53 ± 0.018	0.56 ± 0.03	0.61 ± 0.031	0.52 ± 0.032	0.62 ± 0.043	0.6 ± 0.037	0.58 ± 0.034		
	Shape		0.62 ± 0.029	0.56 ± 0.0085	0.56 ± 0.031	0.63 ± 0.021	0.54 ± 0.015	0.58 ± 0.008	0.54 ± 0.0094	0.64 ± 0.048	0.55 ± 0.0079	0.59 ± 0.012	0.61 ± 0.027	0.52 ± 0.0031	
	WAV	HHH	0.65 ± 0.043	0.55 ± 0.015	0.62 ± 0.0083	0.68 ± 0.062	0.6 ± 0.028	0.55 ± 0.027	0.55 ± 0.026	0.62 ± 0.035	0.51 ± 0.016	0.61 ± 0.052	0.63 ± 0.033	0.56 ± 0.022	
		HHL	0.64 ± 0.054	0.53 ± 0.016	0.67 ± 0.017	0.66 ± 0.083	0.62 ± 0.03	0.59 ± 0.0064	0.55 ± 0.022	0.63 ± 0.017	0.53 ± 0.031	0.62 ± 0.046	0.62 ± 0.083	0.58 ± 0.022	
		HLH	0.66 ± 0.047	0.55 ± 0.035	0.65 ± 0.021	0.68 ± 0.066	0.59 ± 0.024	0.57 ± 0.046	0.61 ± 0.01	0.65 ± 0.025	0.53 ± 0.037	0.64 ± 0.03	0.61 ± 0.054	0.57 ± 0.018	
		HLL	0.65 ± 0.052	0.54 ± 0.015	0.56 ± 0.039	0.67 ± 0.075	0.54 ± 0.0062	0.54 ± 0.027	0.53 ± 0.038	0.61 ± 0.028	0.51 ± 0.019	0.62 ± 0.039	0.59 ± 0.046	0.56 ± 0.021	
		LHH	0.65 ± 0.043	0.54 ± 0.013	0.63 ± 0.017	0.66 ± 0.081	0.59 ± 0.014	0.56 ± 0.0085	0.56 ± 0.03	0.59 ± 0.033	0.52 ± 0.027	0.62 ± 0.037	0.61 ± 0.045	0.55 ± 0.011	
		LHL	0.65 ± 0.049	0.56 ± 0.037	0.7 ± 0.037	0.67 ± 0.075	0.61 ± 0.027	0.54 ± 0.025	0.57 ± 0.045	0.64 ± 0.068	0.52 ± 0.027	0.62 ± 0.043	0.63 ± 0.085	0.62 ± 0.016	
		LLH	0.65 ± 0.047	0.57 ± 0.049	0.63 ± 0.033	0.68 ± 0.063	0.6 ± 0.017	0.58 ± 0.053	0.57 ± 0.025	0.6 ± 0.025	0.52 ± 0.039	0.63 ± 0.037	0.59 ± 0.049	0.54 ± 0.018	
		LLL	0.7 ± 0.039	0.56 ± 0.04	0.56 ± 0.0039	0.69 ± 0.057	0.61 ± 0.038	0.62 ± 0.03	0.59 ± 0.025	0.62 ± 0.051	0.53 ± 0.041	0.66 ± 0.029	0.63 ± 0.063	0.55 ± 0.016	
		CTD	Bin 64		0.62 ± 0.024	0.56 ± 0.039	0.61 ± 0.0083	0.61 ± 0.036	0.61 ± 0.044	0.6 ± 0.055	0.6 ± 0.042	0.62 ± 0.036	0.53 ± 0.036	0.62 ± 0.019	0.71 ± 0.11
	LOG		Sigma: 0.5	0.66 ± 0.033	0.54 ± 0.041	0.66 ± 0.0083	0.64 ± 0.042	0.62 ± 0.041	0.55 ± 0.014	0.6 ± 0.042	0.63 ± 0.027	0.54 ± 0.046	0.65 ± 0.041	0.63 ± 0.072	0.6 ± 0.013
			Sigma: 1	0.66 ± 0.039	0.55 ± 0.035	0.6 ± 0.014	0.66 ± 0.02	0.57 ± 0.023	0.56 ± 0.037	0.63 ± 0.042	0.63 ± 0.035	0.53 ± 0.036	0.65 ± 0.034	0.61 ± 0.04	0.61 ± 0.014
			Sigma: 1.5	0.66 ± 0.049	0.54 ± 0.038	0.57 ± 0.0096	0.65 ± 0.04	0.59 ± 0.014	0.56 ± 0.034	0.63 ± 0.05	0.62 ± 0.037	0.53 ± 0.052	0.66 ± 0.026	0.6 ± 0.053	0.63 ± 0.027
			sigma: 2	0.67 ± 0.025	0.56 ± 0.049	0.62 ± 0.021	0.65 ± 0.03	0.63 ± 0.038	0.56 ± 0.031	0.6 ± 0.019	0.67 ± 0.035	0.52 ± 0.032	0.67 ± 0.034	0.65 ± 0.035	0.63 ± 0.03
			sigma: 2.5	0.65 ± 0.041	0.55 ± 0.041	0.56 ± 0.013	0.64 ± 0.036	0.63 ± 0.05	0.59 ± 0.035	0.58 ± 0.034	0.64 ± 0.033	0.52 ± 0.032	0.65 ± 0.04	0.65 ± 0.057	0.59 ± 0.046
			sigma: 3	0.65 ± 0.047	0.52 ± 0.017	0.56 ± 0.014	0.64 ± 0.041	0.62 ± 0.028	0.56 ± 0.013	0.59 ± 0.05	0.65 ± 0.033	0.52 ± 0.033	0.66 ± 0.028	0.66 ± 0.055	0.62 ± 0.034
sigma: 3.5			0.67 ± 0.024	0.55 ± 0.042	0.57 ± 0.015	0.64 ± 0.037	0.59 ± 0.018	0.6 ± 0.056	0.61 ± 0.041	0.64 ± 0.033	0.52 ± 0.035	0.66 ± 0.029	0.65 ± 0.031	0.61 ± 0.013	
sigma: 4			0.67 ± 0.023	0.53 ± 0.025	0.59 ± 0.04	0.66 ± 0.055	0.6 ± 0.02	0.58 ± 0.014	0.6 ± 0.026	0.64 ± 0.029	0.52 ± 0.028	0.65 ± 0.041	0.65 ± 0.045	0.59 ± 0.014	
sigma: 4.5			0.65 ± 0.041	0.53 ± 0.026	0.64 ± 0.06	0.63 ± 0.051	0.6 ± 0.033	0.57 ± 0	0.57 ± 0.023	0.63 ± 0.024	0.53 ± 0.033	0.63 ± 0.057	0.67 ± 0.018	0.59 ± 0.025	
Sigma: 5	0.66 ± 0.035		0.52 ± 0.014	0.55 ± 0.008	0.64 ± 0.035	0.6 ± 0.031	0.58 ± 0.034	0.59 ± 0.027	0.66 ± 0.057	0.52 ± 0.025	0.64 ± 0.048	0.66 ± 0.035	0.57 ± 0.029		
Shape			0.6 ± 0.025	0.53 ± 0.0041	0.56 ± 0.031	0.59 ± 0.028	0.54 ± 0.0086	0.52 ± 0.0061	0.55 ± 0	0.67 ± 0.084	0.56 ± 0.017	0.58 ± 0.015	0.67 ± 0.032	0.52 ± 0.0068	
WAV	HHH		0.67 ± 0.036	0.61 ± 0.036	0.64 ± 0.0083	0.65 ± 0.032	0.63 ± 0.05	0.65 ± 0.031	0.6 ± 0.025	0.65 ± 0.016	0.54 ± 0.031	0.67 ± 0.02	0.65 ± 0.032	0.63 ± 0.011	
	HHL		0.65 ± 0.046	0.55 ± 0.017	0.67 ± 0	0.71 ± 0.042	0.61 ± 0.04	0.56 ± 0.014	0.57 ± 0.039	0.63 ± 0.039	0.53 ± 0.041	0.66 ± 0.027	0.64 ± 0.025	0.6 ± 0.019	
	HLH		0.66 ± 0.03	0.56 ± 0.046	0.64 ± 0.017	0.63 ± 0.044	0.62 ± 0.058	0.55 ± 0.033	0.62 ± 0.015	0.62 ± 0.033	0.56 ± 0.062	0.65 ± 0.034	0.66 ± 0.039	0.59 ± 0.022	
	HLL		0.65 ± 0.046	0.54 ± 0.04	0.63 ± 0.0004	0.64 ± 0.043	0.58 ± 0.026	0.61 ± 0.041	0.59 ± 0.046	0.66 ± 0.05	0.52 ± 0.028	0.64 ± 0.05	0.61 ± 0.014	0.57 ± 0.017	
	LHH		0.66 ± 0.027	0.6 ± 0.059	0.64 ± 0.0052	0.66 ± 0.021	0.63 ± 0.026	0.58 ± 0.034	0.63 ± 0.037	0.68 ± 0.047	0.53 ± 0.041	0.66 ± 0.029	0.68 ± 0.073	0.61 ± 0.066	
	LHL		0.66 ± 0.039	0.57 ± 0.045	0.67 ± 0.017	0.66 ± 0.02	0.62 ± 0.057	0.56 ± 0.035	0.6 ± 0.035	0.66 ± 0.031	0.52 ± 0.033	0.65 ± 0.037	0.64 ± 0.055	0.63 ± 0.016	
	LLH		0.68 ± 0.045	0.58 ± 0.053	0.63 ± 0.01	0.64 ± 0.055	0.62 ± 0.045	0.63 ± 0.039	0.59 ± 0.05	0.64 ± 0.038	0.56 ± 0.054	0.67 ± 0.03	0.64 ± 0.032	0.63 ± 0.033	
	LLL		0.69 ± 0.031	0.62 ± 0.084	0.66 ± 0.053	0.65 ± 0.024	0.65 ± 0.09	0.64 ± 0.053	0.67 ± 0.11	0.71 ± 0.056	0.53 ± 0.045	0.68 ± 0.027	0.72 ± 0.074	0.7 ± 0.069	
	PET		Bin 64		0.74 ± 0.018	0.53 ± 0.026	0.75 ± 0.0096	0.68 ± 0.032	0.59 ± 0.02	0.57 ± 0.039	0.59 ± 0.059	0.66 ± 0.1	0.54 ± 0.062	0.71 ± 0.017	0.64 ± 0.11
LOG			Sigma: 0.5	0.66 ± 0.021	0.55 ± 0.044	0.7 ± 0.019	0.64 ± 0.023	0.63 ± 0.073	0.63 ± 0.018	0.59 ± 0.034	0.63 ± 0.08	0.52 ± 0.039	0.65 ± 0.057	0.6 ± 0.092	0.61 ± 0.023
			Sigma: 1	0.64 ± 0.012	0.56 ± 0.05	0.74 ± 0.03	0.62 ± 0.037	0.63 ± 0.073	0.63 ± 0.022	0.63 ± 0.079	0.6 ± 0.051	0.55 ± 0.048	0.64 ± 0.061	0.62 ± 0.099	0.66 ± 0.058
			Sigma: 1.5	0.67 ± 0.019	0.52 ± 0.029	0.74 ± 0.019	0.65 ± 0.025	0.63 ± 0.046	0.64 ± 0.028	0.57 ± 0.048	0.58 ± 0.056	0.55 ± 0.05	0.65 ± 0.045	0.63 ± 0.1	0.56 ± 0.034
			sigma: 2	0.64 ± 0.02	0.54 ± 0.03	0.73 ± 0.026	0.64 ± 0.024	0.6 ± 0.05	0.63 ± 0.025	0.56 ± 0.062	0.61 ± 0.048	0.52 ± 0.016	0.66 ± 0.039	0.67 ± 0.13	0.59 ± 0.022
			sigma: 2.5	0.71 ± 0.036	0.59 ± 0.072	0.67 ± 0.053	0.7 ± 0.05	0.59 ± 0.018	0.66 ± 0.015	0.63 ± 0.085	0.63 ± 0.1	0.54 ± 0.045	0.69 ± 0.046	0.66 ± 0.12	0.63 ± 0.069
			sigma: 3	0.62 ± 0.033	0.58 ± 0.072	0.71 ± 0.031	0.65 ± 0.018	0.58 ± 0.022	0.6 ± 0.047	0.54 ± 0.032	0.59 ± 0.084	0.53 ± 0.035	0.65 ± 0.047	0.64 ± 0.11	0.57 ± 0.023
		sigma: 3.5	0.64 ± 0.023	0.54 ± 0.074	0.71 ± 0.012	0.65 ± 0.037	0.58 ± 0.019	0.6 ± 0.06	0.56 ± 0.048	0.6 ± 0.024	0.54 ± 0.037	0.65 ± 0.056	0.69 ± 0.15	0.62 ± 0.033	
		sigma: 4	0.63 ± 0.026	0.57 ± 0.068	0.73 ± 0.012	0.66 ± 0.013	0.61 ± 0.034	0.57 ± 0.022	0.55 ± 0.031	0.65 ± 0.097	0.53 ± 0.031	0.64 ± 0.061	0.67 ± 0.13	0.56 ± 0.017	
		sigma: 4.5	0.63 ± 0.02	0.54 ± 0.033	0.7 ± 0.0084	0.63 ± 0.032	0.58 ± 0.03	0.55 ± 0.0004	0.54 ± 0.043	0.57 ± 0.049	0.56 ± 0.064	0.65 ± 0.048	0.63 ± 0.11	0.55 ± 0.018	
Sigma: 5		0.63 ± 0.035	0.6 ± 0.083	0.7 ± 0.016	0.62 ± 0.037	0.6 ± 0.012	0.55 ± 0.0088	0.55 ± 0.036	0.62 ± 0.069	0.53 ± 0.033	0.66 ± 0.046	0.6 ± 0.084	0.63 ± 0.0086		
Shape			0.58 ± 0.024	0.51 ± 0.0047	0.6 ± 0.05	0.64 ± 0.019	0.55 ± 0.019	0.54 ± 0.0046	0.52 ± 0	0.7 ± 0.066	0.56 ± 0.012	0.58 ± 0.021	0.68 ± 0.046	0.52 ± 0.0037	
WAV		HHH	0.67 ± 0.026	0.69 ± 0.032	0.65 ± 0.0078	0.68 ± 0.021	0.66 ± 0.02	0.72 ± 0.016	0.6 ± 0.043	0.61 ± 0.063	0.56 ± 0.075	0.73 ± 0.013	0.65 ± 0.12	0.64 ± 0.038	
		HHL	0.65 ± 0.014	0.62 ± 0.081	0.69 ± 0.016	0.63 ± 0.038	0.66 ± 0.014	0.57 ± 0.031	0.59 ± 0.056	0.63 ± 0.066	0.52 ± 0.024	0.7 ± 0.0094	0.64 ± 0.11	0.65 ± 0.014	
		HLH	0.63 ± 0.022	0.62 ± 0.045	0.69 ± 0.019	0.62 ± 0.04	0.66 ± 0.014	0.62 ± 0	0.56 ± 0.043	0.57 ± 0.041	0.53 ± 0.03	0.67 ± 0.027	0.62 ± 0.092	0.57 ± 0.018	
		HLL	0.63 ± 0.03	0.54 ± 0.02											

- **KRAS**

Fig 3 shows a heatmap of KRAS mutation status prediction different combination of feature selection, classifier and image sets. The model performance has a wide range from 0.5 to 0.83, which is as same as EGFR models. According to this figure, the combination methods with highest predictive performances are: LOG preprocessed of CT image set with sigma 3.5 with SM feature selector and SGD classifier (CT_LOG_3.5S+SM+SGD, AUC: 0.83), CT_LOG_5 S image set with VT_SM feature selector SGD classifier (CT_LOG_5 S+ VT_SM+SGD, AUC: 0.82). These models are followed by CTD_LOG_4 S image set with VT feature selector and SGD classifier (CTD_LOG_4S+VT+SGD, AUC: 0.80), PET_LOG_0 S image set with SM feature selector and SGD classifier (PET_LOG_0.5S+SM+SGD, AUC: 0.80) and PET image set which preprocessed by HHH of wavelet with SKB feature selector and BNB classifier (PET_W_HHH+SKB+BNB, AUC: 0.80).

Supplemental Tables 8 and 9 show KRAS mutation status prediction results based on feature selection methods and image sets. In these results, feature selection performance has a range from 0.5 to 0.83, and the combination of SM feature selection with PET_LOG_4.5S image set had the highest performance (AUC: 0.66 ± 0.093 & 0.51 - 0.79), followed by VT_SM feature selection with CT_LOG_3.5S image set (AUC: 0.66 ± 0.079 & 0.55 - 0.8), SM feature selection with CT_LOG_3.5S image set (AUC: 0.65 ± 0.11 & 0.5 - 0.83) and SM feature selection with PET_BIN64 image set (AUC: 0.64 ± 0.058 & 0.51 - 0.71). Supplemental figure 12 delicate the box plot of KRAS mutation status prediction based on feature selection methods.

KRAS mutation status prediction results based on classifier and image sets is presented in supplemental Tables 10 and 11 (mean \pm SD & Min-Max). Here, the classifier performance has a range from 0.5 to 0.83 and combination of SVM classifier with CT_LOG_3.5S image set had the highest performance (AUC: 0.79 ± 0.021 & 0.76-0.8), followed by BNB classifier PET_W_HHH image set (AUC: 0.75 ± 0.1 & 0.6 - 0.81) and SGD classifier with CTD_LOG_4 S image set (AUC: 0.71 ± 0.08 & 0.56 - 0.81). Supplemental figure 13 delicate the box plot of KRAS mutation status prediction based on feature selection methods.

KRAS AUC in Independent Validation Set




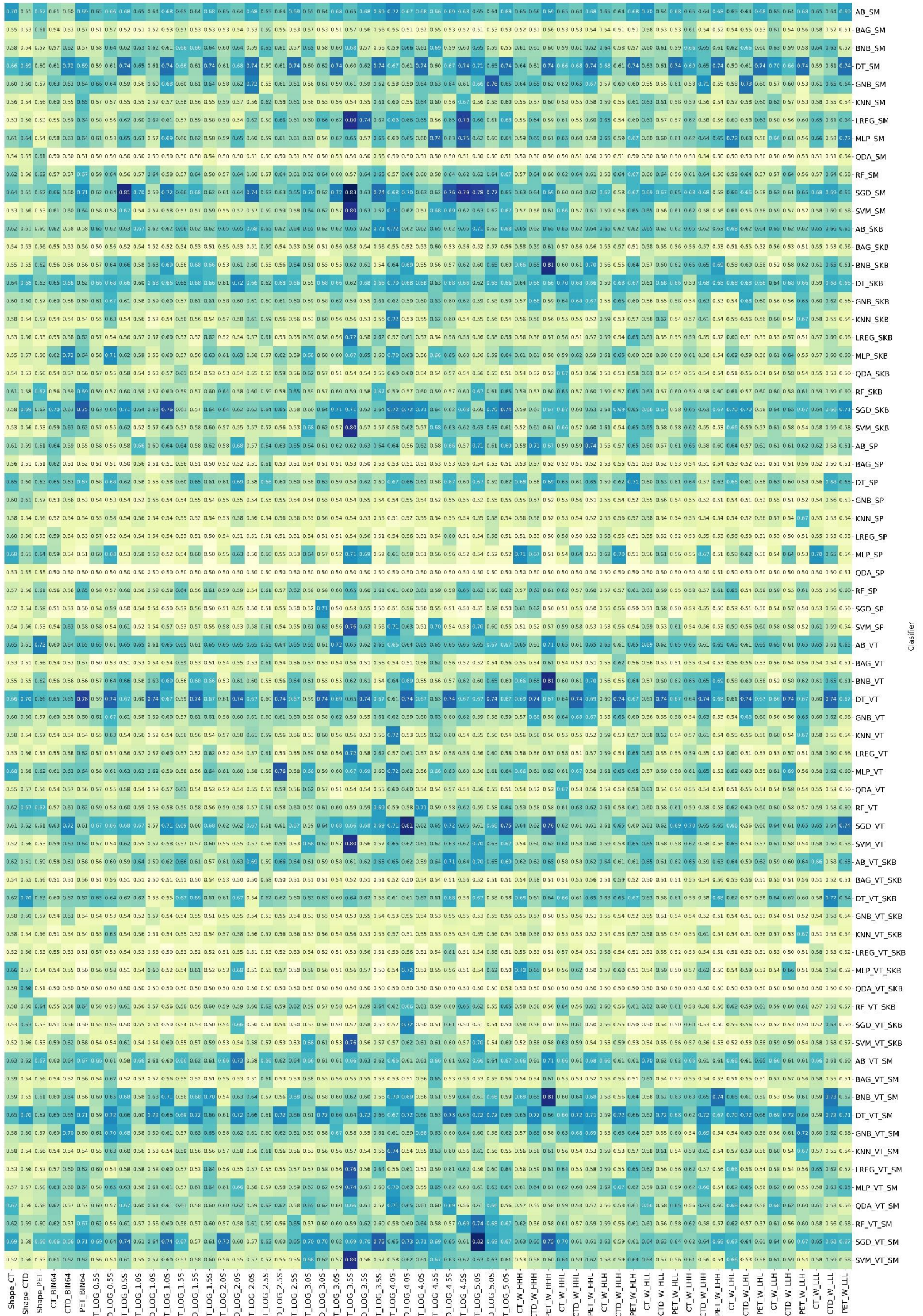


Figure 3. Heatmap depicting the predictive performance (AUC) of feature selection-classification (rows) and image sets (columns) in prediction of KRAS mutation status of independent validation sets

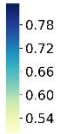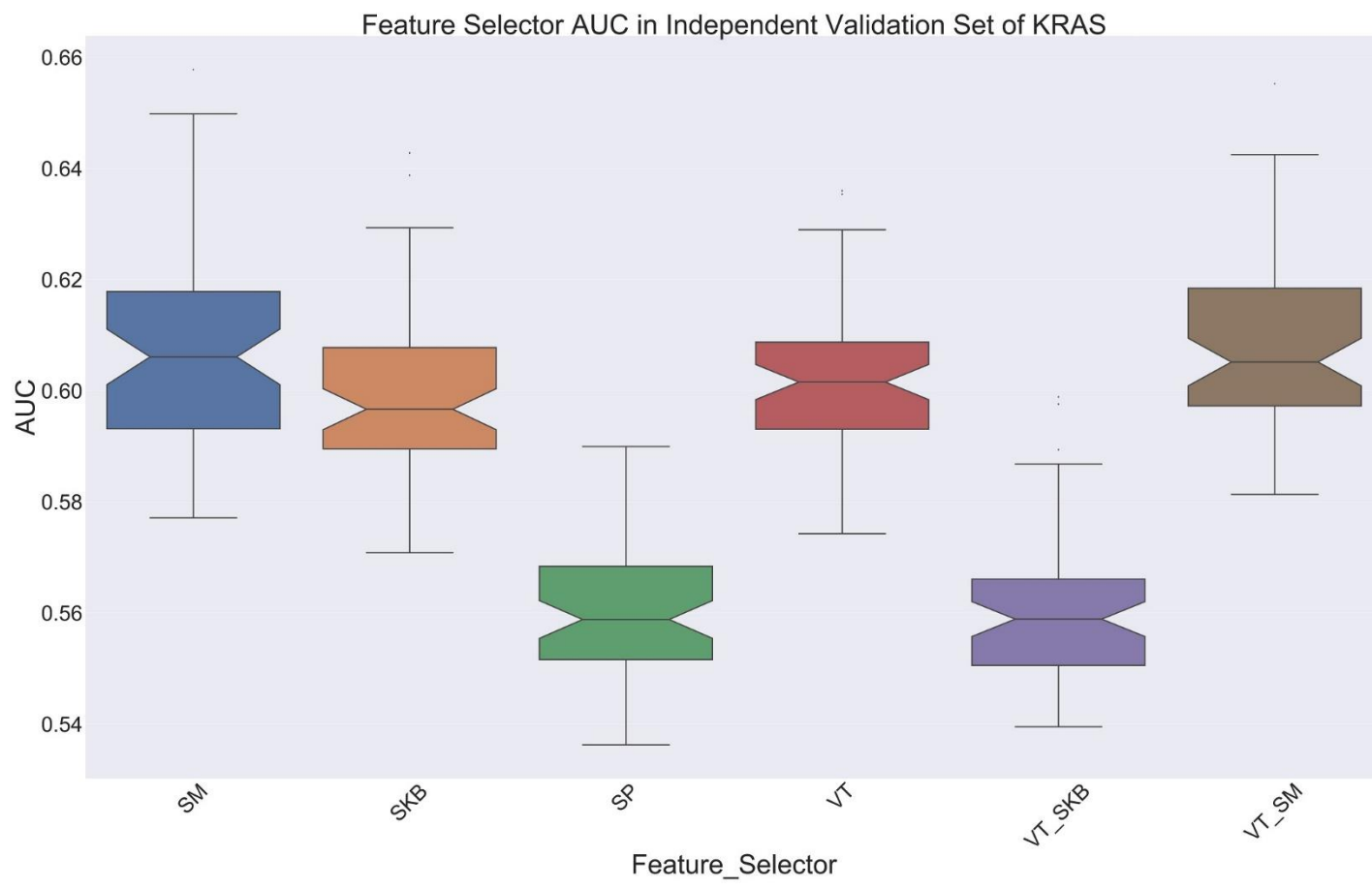

Supplemental Fig 12. AUC box plot of different feature selector in independent validation set of KRAS

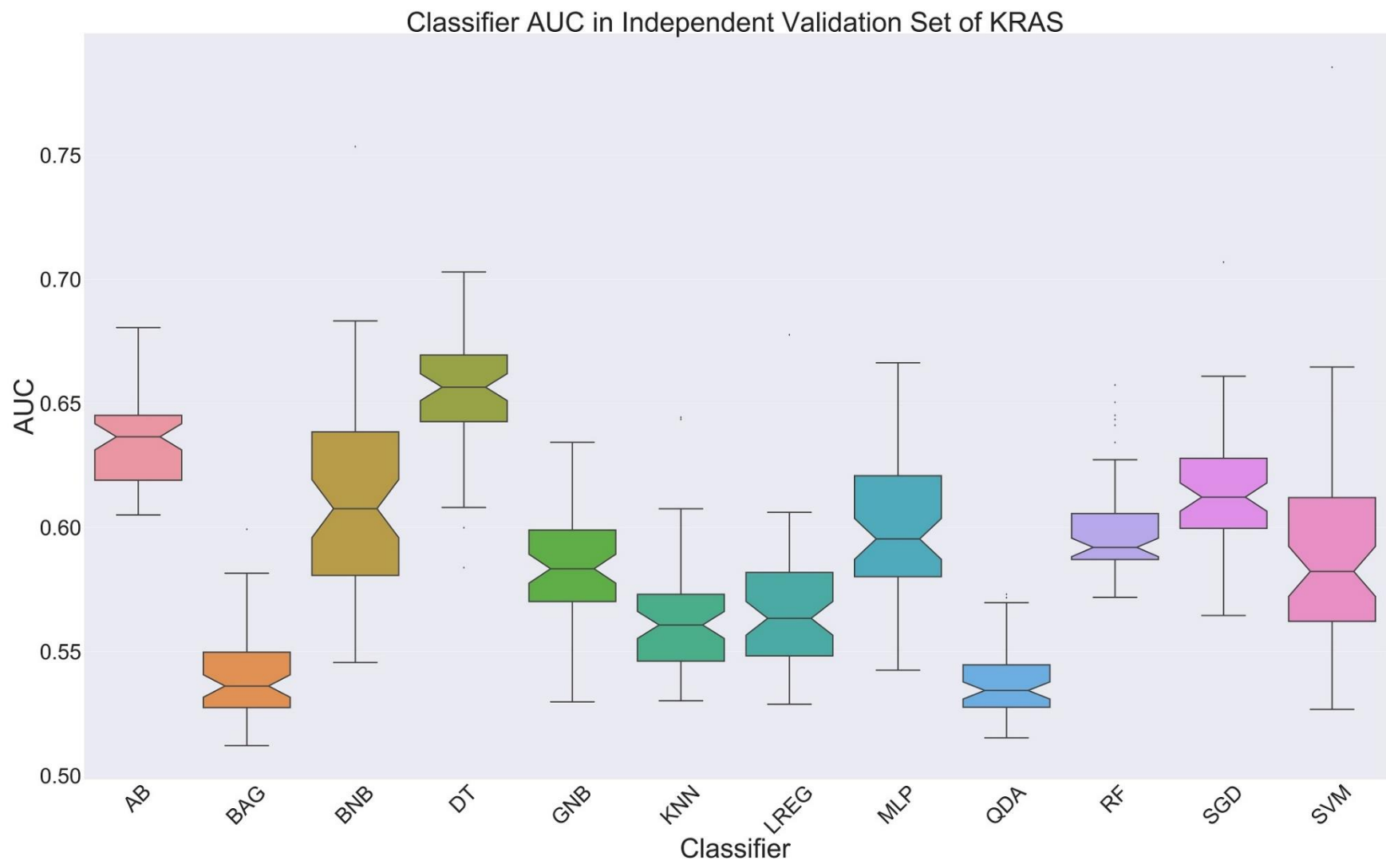

Supplemental Fig 13. AUC box plot of different classifier in independent validation set of KRAS

Supplemental Table 8. KRAS Feature Selector Min-Max different modality, preprocessing and setting

Modality	Preprocessed	Setting	SM	SKB	SP	VT	VT_SKB	VT_SM
CT	Bin 64		0.5 - 0.66	0.54 - 0.7	0.5 - 0.65	0.54 - 0.65	0.5 - 0.6	0.54 - 0.66
	LOG	Sigma: 0.5	0.5 - 0.66	0.5 - 0.66	0.5 - 0.6	0.5 - 0.67	0.5 - 0.65	0.54 - 0.69
		Sigma: 1	0.5 - 0.7	0.52 - 0.67	0.5 - 0.66	0.53 - 0.67	0.5 - 0.64	0.53 - 0.66
		Sigma: 1.5	0.5 - 0.66	0.54 - 0.66	0.5 - 0.64	0.56 - 0.69	0.5 - 0.67	0.55 - 0.69
		sigma: 2	0.5 - 0.65	0.52 - 0.65	0.5 - 0.61	0.52 - 0.65	0.5 - 0.61	0.54 - 0.73
		sigma: 2.5	0.5 - 0.65	0.55 - 0.65	0.5 - 0.66	0.55 - 0.65	0.5 - 0.62	0.55 - 0.66
		sigma: 3	0.5 - 0.7	0.53 - 0.68	0.5 - 0.64	0.53 - 0.68	0.5 - 0.68	0.53 - 0.7
		sigma: 3.5	0.5 - 0.83	0.54 - 0.8	0.5 - 0.76	0.54 - 0.8	0.5 - 0.76	0.55 - 0.8
		sigma: 4	0.5 - 0.72	0.55 - 0.72	0.5 - 0.71	0.54 - 0.72	0.5 - 0.65	0.53 - 0.74
		sigma: 4.5	0.5 - 0.74	0.54 - 0.68	0.5 - 0.7	0.54 - 0.66	0.5 - 0.64	0.58 - 0.69
	Sigma: 5	0.5 - 0.78	0.52 - 0.71	0.5 - 0.71	0.52 - 0.7	0.5 - 0.7	0.52 - 0.82	
	Shape		0.53 - 0.7	0.53 - 0.64	0.52 - 0.68	0.52 - 0.68	0.52 - 0.66	0.52 - 0.69
	WAV	HHH	0.5 - 0.65	0.52 - 0.66	0.5 - 0.71	0.54 - 0.69	0.5 - 0.7	0.53 - 0.68
		HHL	0.5 - 0.66	0.55 - 0.7	0.5 - 0.65	0.54 - 0.67	0.5 - 0.66	0.55 - 0.7
		HLH	0.5 - 0.68	0.53 - 0.65	0.5 - 0.62	0.53 - 0.65	0.5 - 0.64	0.55 - 0.66
		HLL	0.5 - 0.7	0.56 - 0.66	0.5 - 0.61	0.56 - 0.69	0.5 - 0.65	0.57 - 0.7
		LHH	0.5 - 0.69	0.54 - 0.65	0.5 - 0.65	0.54 - 0.7	0.5 - 0.65	0.54 - 0.66
		LHL	0.5 - 0.72	0.51 - 0.7	0.5 - 0.66	0.51 - 0.66	0.5 - 0.64	0.51 - 0.7
		LLH	0.5 - 0.7	0.52 - 0.65	0.5 - 0.63	0.52 - 0.66	0.5 - 0.64	0.53 - 0.69
		LLL	0.51 - 0.68	0.55 - 0.65	0.5 - 0.7	0.54 - 0.65	0.5 - 0.66	0.57 - 0.66
CTD	Bin 64		0.5 - 0.72	0.53 - 0.72	0.5 - 0.63	0.53 - 0.72	0.5 - 0.62	0.52 - 0.7
	LOG	Sigma: 0.5	0.5 - 0.66	0.54 - 0.71	0.5 - 0.68	0.52 - 0.74	0.5 - 0.64	0.54 - 0.72
		Sigma: 1	0.5 - 0.64	0.52 - 0.68	0.5 - 0.6	0.52 - 0.74	0.5 - 0.6	0.52 - 0.72
		Sigma: 1.5	0.5 - 0.68	0.52 - 0.68	0.5 - 0.65	0.52 - 0.74	0.5 - 0.69	0.52 - 0.72
		sigma: 2	0.5 - 0.68	0.53 - 0.72	0.5 - 0.69	0.54 - 0.74	0.5 - 0.68	0.53 - 0.73
		sigma: 2.5	0.5 - 0.66	0.53 - 0.68	0.5 - 0.63	0.53 - 0.76	0.5 - 0.66	0.53 - 0.72
		sigma: 3	0.5 - 0.66	0.51 - 0.68	0.5 - 0.71	0.55 - 0.74	0.5 - 0.63	0.55 - 0.72
		sigma: 3.5	0.5 - 0.74	0.53 - 0.68	0.5 - 0.69	0.51 - 0.74	0.5 - 0.62	0.54 - 0.72
		sigma: 4	0.5 - 0.7	0.52 - 0.72	0.5 - 0.66	0.51 - 0.81	0.5 - 0.72	0.51 - 0.73
		sigma: 4.5	0.5 - 0.76	0.52 - 0.68	0.5 - 0.67	0.56 - 0.74	0.5 - 0.71	0.53 - 0.73
	Sigma: 5	0.5 - 0.77	0.55 - 0.7	0.5 - 0.61	0.54 - 0.74	0.5 - 0.65	0.55 - 0.72	
	Shape		0.53 - 0.69	0.53 - 0.69	0.51 - 0.61	0.51 - 0.7	0.54 - 0.7	0.54 - 0.7
	WAV	HHH	0.5 - 0.66	0.52 - 0.68	0.5 - 0.71	0.52 - 0.74	0.5 - 0.65	0.54 - 0.72
		HHL	0.51 - 0.68	0.51 - 0.68	0.5 - 0.64	0.51 - 0.74	0.5 - 0.62	0.53 - 0.72
		HLH	0.51 - 0.64	0.53 - 0.69	0.5 - 0.7	0.53 - 0.74	0.5 - 0.65	0.53 - 0.72
		HLL	0.5 - 0.67	0.54 - 0.68	0.5 - 0.63	0.53 - 0.74	0.5 - 0.61	0.54 - 0.72
		LHH	0.54 - 0.71	0.54 - 0.68	0.5 - 0.67	0.54 - 0.74	0.5 - 0.62	0.54 - 0.72
		LHL	0.5 - 0.73	0.51 - 0.7	0.5 - 0.62	0.51 - 0.74	0.5 - 0.59	0.53 - 0.72
		LLH	0.5 - 0.66	0.53 - 0.68	0.5 - 0.64	0.54 - 0.74	0.5 - 0.66	0.54 - 0.72
		LLL	0.51 - 0.69	0.53 - 0.68	0.5 - 0.68	0.53 - 0.74	0.5 - 0.72	0.55 - 0.73
PET	Bin 64		0.51 - 0.71	0.54 - 0.75	0.5 - 0.67	0.54 - 0.78	0.5 - 0.64	0.54 - 0.71
	LOG	Sigma: 0.5	0.5 - 0.81	0.52 - 0.71	0.5 - 0.62	0.51 - 0.68	0.5 - 0.62	0.52 - 0.74
		Sigma: 1	0.5 - 0.74	0.52 - 0.76	0.5 - 0.64	0.54 - 0.71	0.5 - 0.62	0.54 - 0.74
		Sigma: 1.5	0.53 - 0.74	0.51 - 0.66	0.5 - 0.62	0.51 - 0.68	0.5 - 0.61	0.51 - 0.7
		sigma: 2	0.51 - 0.74	0.51 - 0.68	0.5 - 0.58	0.53 - 0.67	0.5 - 0.69	0.51 - 0.66
		sigma: 2.5	0.51 - 0.74	0.52 - 0.66	0.5 - 0.65	0.53 - 0.67	0.5 - 0.64	0.52 - 0.68
		sigma: 3	0.51 - 0.74	0.51 - 0.71	0.5 - 0.62	0.51 - 0.72	0.5 - 0.6	0.56 - 0.67
		sigma: 3.5	0.56 - 0.74	0.54 - 0.71	0.5 - 0.64	0.54 - 0.69	0.5 - 0.65	0.53 - 0.75
		sigma: 4	0.51 - 0.74	0.53 - 0.71	0.5 - 0.62	0.54 - 0.71	0.5 - 0.62	0.51 - 0.71
		sigma: 4.5	0.51 - 0.79	0.54 - 0.68	0.5 - 0.65	0.52 - 0.67	0.5 - 0.65	0.56 - 0.69
	Sigma: 5	0.51 - 0.74	0.51 - 0.74	0.5 - 0.69	0.51 - 0.75	0.5 - 0.69	0.55 - 0.67	
	Shape		0.53 - 0.67	0.53 - 0.67	0.51 - 0.64	0.53 - 0.72	0.51 - 0.64	0.53 - 0.67
	WAV	HHH	0.5 - 0.74	0.53 - 0.81	0.5 - 0.69	0.53 - 0.81	0.5 - 0.65	0.58 - 0.81
		HHL	0.5 - 0.74	0.52 - 0.7	0.5 - 0.74	0.51 - 0.7	0.5 - 0.65	0.52 - 0.71
		HLH	0.5 - 0.74	0.51 - 0.67	0.5 - 0.71	0.56 - 0.67	0.5 - 0.67	0.51 - 0.66
		HLL	0.5 - 0.74	0.54 - 0.66	0.5 - 0.61	0.51 - 0.69	0.5 - 0.62	0.52 - 0.68
		LHH	0.5 - 0.74	0.52 - 0.69	0.5 - 0.63	0.52 - 0.69	0.5 - 0.68	0.53 - 0.74
		LHL	0.5 - 0.74	0.52 - 0.66	0.5 - 0.61	0.53 - 0.67	0.5 - 0.62	0.51 - 0.66
		LLH	0.52 - 0.74	0.51 - 0.67	0.5 - 0.67	0.51 - 0.67	0.5 - 0.67	0.54 - 0.72
	LLL	0.51 - 0.74	0.5 - 0.71	0.5 - 0.65	0.5 - 0.74	0.5 - 0.65	0.51 - 0.71	

Supplemental Table 9. KRAS Feature Selector Mean+Sd different modality, preprocessing and setting

Modality	Preprocessed	Setting	SM	SKB	SP	VT	VT_SKB	VT_SM
CT	Bin 64		0.59 ± 0.043	0.59 ± 0.049	0.57 ± 0.053	0.58 ± 0.036	0.55 ± 0.035	0.6 ± 0.039
	LOG	Sigma: 0.5	0.59 ± 0.046	0.59 ± 0.045	0.56 ± 0.033	0.59 ± 0.044	0.55 ± 0.04	0.6 ± 0.044
		Sigma: 1	0.59 ± 0.062	0.58 ± 0.042	0.56 ± 0.044	0.59 ± 0.043	0.55 ± 0.043	0.58 ± 0.034
		Sigma: 1.5	0.6 ± 0.05	0.59 ± 0.036	0.56 ± 0.047	0.6 ± 0.038	0.56 ± 0.055	0.6 ± 0.046
		sigma: 2	0.58 ± 0.04	0.58 ± 0.045	0.55 ± 0.036	0.57 ± 0.041	0.55 ± 0.038	0.6 ± 0.054
		sigma: 2.5	0.59 ± 0.043	0.59 ± 0.032	0.58 ± 0.055	0.59 ± 0.029	0.56 ± 0.041	0.6 ± 0.031
		sigma: 3	0.61 ± 0.053	0.61 ± 0.046	0.57 ± 0.047	0.61 ± 0.045	0.57 ± 0.056	0.62 ± 0.051
		sigma: 3.5	0.65 ± 0.11	0.64 ± 0.075	0.59 ± 0.084	0.64 ± 0.074	0.57 ± 0.075	0.66 ± 0.079
		sigma: 4	0.62 ± 0.069	0.64 ± 0.064	0.58 ± 0.073	0.64 ± 0.06	0.56 ± 0.049	0.64 ± 0.067
		sigma: 4.5	0.62 ± 0.068	0.61 ± 0.048	0.56 ± 0.052	0.61 ± 0.043	0.56 ± 0.046	0.63 ± 0.032
	Sigma: 5	0.63 ± 0.079	0.6 ± 0.054	0.58 ± 0.078	0.6 ± 0.052	0.58 ± 0.076	0.64 ± 0.082	
	Shape		0.59 ± 0.056	0.57 ± 0.039	0.59 ± 0.05	0.59 ± 0.056	0.58 ± 0.046	0.6 ± 0.052
	WAV	HHH	0.58 ± 0.048	0.59 ± 0.044	0.57 ± 0.07	0.6 ± 0.056	0.58 ± 0.063	0.6 ± 0.05
		HHL	0.61 ± 0.051	0.63 ± 0.047	0.56 ± 0.048	0.61 ± 0.039	0.57 ± 0.051	0.62 ± 0.043
		HLH	0.61 ± 0.055	0.59 ± 0.033	0.56 ± 0.033	0.59 ± 0.033	0.57 ± 0.045	0.59 ± 0.033
		HLL	0.61 ± 0.056	0.61 ± 0.035	0.56 ± 0.039	0.6 ± 0.04	0.56 ± 0.055	0.62 ± 0.038
		LHH	0.62 ± 0.05	0.6 ± 0.038	0.56 ± 0.046	0.6 ± 0.048	0.56 ± 0.04	0.6 ± 0.037
		LHL	0.6 ± 0.057	0.6 ± 0.065	0.57 ± 0.055	0.6 ± 0.047	0.56 ± 0.049	0.63 ± 0.059
		LLH	0.61 ± 0.052	0.59 ± 0.043	0.56 ± 0.041	0.59 ± 0.048	0.56 ± 0.041	0.61 ± 0.048
		LLL	0.61 ± 0.049	0.59 ± 0.031	0.57 ± 0.058	0.59 ± 0.033	0.56 ± 0.047	0.6 ± 0.031
CTD	Bin 64		0.59 ± 0.055	0.6 ± 0.057	0.56 ± 0.045	0.6 ± 0.054	0.57 ± 0.045	0.61 ± 0.053
	LOG	Sigma: 0.5	0.6 ± 0.046	0.61 ± 0.057	0.57 ± 0.06	0.61 ± 0.065	0.57 ± 0.043	0.62 ± 0.058
		Sigma: 1	0.59 ± 0.043	0.58 ± 0.05	0.55 ± 0.031	0.59 ± 0.06	0.55 ± 0.035	0.6 ± 0.052
		Sigma: 1.5	0.6 ± 0.052	0.58 ± 0.052	0.56 ± 0.045	0.59 ± 0.063	0.56 ± 0.057	0.6 ± 0.057
		sigma: 2	0.6 ± 0.055	0.6 ± 0.054	0.58 ± 0.068	0.6 ± 0.057	0.59 ± 0.065	0.61 ± 0.063
		sigma: 2.5	0.6 ± 0.042	0.59 ± 0.045	0.56 ± 0.038	0.61 ± 0.073	0.56 ± 0.049	0.59 ± 0.051
		sigma: 3	0.59 ± 0.049	0.59 ± 0.043	0.58 ± 0.065	0.6 ± 0.052	0.56 ± 0.04	0.61 ± 0.052
		sigma: 3.5	0.61 ± 0.066	0.59 ± 0.046	0.58 ± 0.06	0.6 ± 0.071	0.56 ± 0.039	0.6 ± 0.056
		sigma: 4	0.61 ± 0.063	0.61 ± 0.061	0.57 ± 0.048	0.63 ± 0.084	0.6 ± 0.08	0.62 ± 0.065
		sigma: 4.5	0.63 ± 0.074	0.6 ± 0.041	0.57 ± 0.056	0.62 ± 0.058	0.59 ± 0.064	0.62 ± 0.053
	Sigma: 5	0.62 ± 0.078	0.61 ± 0.047	0.57 ± 0.038	0.62 ± 0.06	0.56 ± 0.043	0.63 ± 0.052	
	Shape		0.58 ± 0.048	0.58 ± 0.052	0.57 ± 0.032	0.59 ± 0.055	0.6 ± 0.049	0.58 ± 0.045
	WAV	HHH	0.6 ± 0.054	0.61 ± 0.047	0.59 ± 0.066	0.61 ± 0.062	0.58 ± 0.041	0.61 ± 0.052
		HHL	0.59 ± 0.051	0.59 ± 0.053	0.56 ± 0.046	0.6 ± 0.07	0.57 ± 0.036	0.6 ± 0.056
		HLH	0.58 ± 0.042	0.59 ± 0.061	0.57 ± 0.052	0.6 ± 0.064	0.57 ± 0.044	0.59 ± 0.054
		HLL	0.59 ± 0.049	0.59 ± 0.046	0.56 ± 0.041	0.59 ± 0.057	0.55 ± 0.037	0.6 ± 0.048
		LHH	0.62 ± 0.054	0.6 ± 0.042	0.56 ± 0.045	0.61 ± 0.058	0.56 ± 0.039	0.62 ± 0.054
		LHL	0.61 ± 0.062	0.6 ± 0.059	0.55 ± 0.04	0.59 ± 0.063	0.55 ± 0.035	0.6 ± 0.05
		LLH	0.59 ± 0.045	0.59 ± 0.045	0.56 ± 0.047	0.6 ± 0.06	0.56 ± 0.055	0.59 ± 0.048
		LLL	0.6 ± 0.054	0.6 ± 0.05	0.57 ± 0.056	0.61 ± 0.057	0.58 ± 0.06	0.62 ± 0.059
PET	Bin 64		0.64 ± 0.058	0.61 ± 0.062	0.55 ± 0.058	0.62 ± 0.063	0.55 ± 0.05	0.62 ± 0.059
	LOG	Sigma: 0.5	0.63 ± 0.09	0.61 ± 0.054	0.55 ± 0.045	0.61 ± 0.053	0.55 ± 0.046	0.63 ± 0.059
		Sigma: 1	0.62 ± 0.077	0.61 ± 0.067	0.55 ± 0.042	0.61 ± 0.057	0.54 ± 0.039	0.62 ± 0.056
		Sigma: 1.5	0.61 ± 0.061	0.6 ± 0.05	0.54 ± 0.049	0.61 ± 0.058	0.55 ± 0.041	0.61 ± 0.049
		sigma: 2	0.62 ± 0.08	0.6 ± 0.046	0.54 ± 0.033	0.6 ± 0.046	0.55 ± 0.059	0.59 ± 0.038
		sigma: 2.5	0.61 ± 0.063	0.59 ± 0.048	0.55 ± 0.052	0.6 ± 0.049	0.55 ± 0.051	0.61 ± 0.05
		sigma: 3	0.61 ± 0.071	0.59 ± 0.055	0.54 ± 0.042	0.6 ± 0.067	0.54 ± 0.035	0.6 ± 0.037
		sigma: 3.5	0.64 ± 0.063	0.61 ± 0.054	0.55 ± 0.047	0.61 ± 0.057	0.55 ± 0.056	0.6 ± 0.059
		sigma: 4	0.62 ± 0.064	0.59 ± 0.057	0.54 ± 0.048	0.6 ± 0.055	0.55 ± 0.047	0.6 ± 0.059
		sigma: 4.5	0.66 ± 0.093	0.61 ± 0.046	0.55 ± 0.046	0.6 ± 0.049	0.55 ± 0.052	0.61 ± 0.043
	Sigma: 5	0.63 ± 0.07	0.61 ± 0.063	0.56 ± 0.063	0.62 ± 0.066	0.56 ± 0.065	0.62 ± 0.046	
	Shape		0.59 ± 0.043	0.58 ± 0.043	0.57 ± 0.044	0.6 ± 0.06	0.56 ± 0.041	0.59 ± 0.045
	WAV	HHH	0.62 ± 0.065	0.62 ± 0.072	0.55 ± 0.071	0.64 ± 0.085	0.55 ± 0.056	0.64 ± 0.074
		HHL	0.61 ± 0.062	0.61 ± 0.052	0.55 ± 0.076	0.61 ± 0.062	0.55 ± 0.056	0.61 ± 0.064
		HLH	0.62 ± 0.07	0.62 ± 0.046	0.56 ± 0.07	0.62 ± 0.039	0.55 ± 0.056	0.61 ± 0.04
		HLL	0.6 ± 0.065	0.59 ± 0.037	0.55 ± 0.042	0.59 ± 0.056	0.55 ± 0.048	0.6 ± 0.048
		LHH	0.59 ± 0.07	0.58 ± 0.068	0.55 ± 0.044	0.58 ± 0.067	0.55 ± 0.059	0.59 ± 0.068
		LHL	0.6 ± 0.064	0.58 ± 0.042	0.54 ± 0.044	0.58 ± 0.046	0.55 ± 0.043	0.58 ± 0.042
		LLH	0.59 ± 0.069	0.6 ± 0.063	0.56 ± 0.057	0.6 ± 0.06	0.56 ± 0.06	0.61 ± 0.057
		LLL	0.62 ± 0.075	0.59 ± 0.059	0.55 ± 0.048	0.59 ± 0.068	0.55 ± 0.053	0.6 ± 0.056

Supplemental Table 10. KRAS Classifier Min-Max different modality, preprocessing and setting

Modality	Preprocessed	Setting	AB	BAG	BNB	DT	GNB	KNN	LREG	MLP	QDA	RF	SGD	SVM	
CT	Bin 64		0.58 - 0.64	0.51 - 0.62	0.56 - 0.6	0.6 - 0.65	0.53 - 0.63	0.51 - 0.6	0.55 - 0.59	0.54 - 0.63	0.5 - 0.62	0.55 - 0.62	0.51 - 0.7	0.54 - 0.61	
	LOG	Sigma: 0.5	0.56 - 0.66	0.5 - 0.54	0.57 - 0.6	0.58 - 0.66	0.54 - 0.66	0.55 - 0.57	0.57 - 0.6	0.56 - 0.61	0.5 - 0.6	0.58 - 0.62	0.54 - 0.69	0.54 - 0.58	
		Sigma: 1	0.64 - 0.67	0.51 - 0.56	0.58 - 0.63	0.58 - 0.65	0.52 - 0.58	0.55 - 0.56	0.53 - 0.6	0.54 - 0.63	0.5 - 0.6	0.56 - 0.59	0.54 - 0.7	0.52 - 0.55	
		Sigma: 1.5	0.64 - 0.66	0.51 - 0.59	0.56 - 0.66	0.59 - 0.69	0.54 - 0.6	0.55 - 0.59	0.53 - 0.57	0.54 - 0.6	0.5 - 0.61	0.57 - 0.64	0.53 - 0.69	0.53 - 0.57	
		sigma: 2	0.58 - 0.66	0.52 - 0.55	0.53 - 0.6	0.61 - 0.61	0.55 - 0.58	0.53 - 0.59	0.51 - 0.58	0.53 - 0.61	0.5 - 0.6	0.57 - 0.64	0.51 - 0.73	0.53 - 0.58	
		sigma: 2.5	0.59 - 0.66	0.58 - 0.61	0.55 - 0.65	0.59 - 0.66	0.53 - 0.6	0.54 - 0.62	0.51 - 0.61	0.55 - 0.6	0.5 - 0.62	0.57 - 0.62	0.51 - 0.64	0.56 - 0.61	
		sigma: 3	0.61 - 0.66	0.51 - 0.58	0.61 - 0.65	0.58 - 0.63	0.53 - 0.59	0.53 - 0.56	0.54 - 0.6	0.58 - 0.68	0.5 - 0.65	0.57 - 0.62	0.52 - 0.7	0.61 - 0.68	
		sigma: 3.5	0.61 - 0.66	0.52 - 0.58	0.62 - 0.68	0.58 - 0.65	0.53 - 0.61	0.54 - 0.57	0.53 - 0.8	0.56 - 0.74	0.5 - 0.66	0.54 - 0.65	0.52 - 0.83	0.76 - 0.8	
		sigma: 4	0.64 - 0.72	0.52 - 0.55	0.59 - 0.7	0.61 - 0.7	0.53 - 0.62	0.51 - 0.74	0.53 - 0.68	0.54 - 0.72	0.5 - 0.71	0.57 - 0.62	0.51 - 0.72	0.57 - 0.71	
		sigma: 4.5	0.58 - 0.66	0.53 - 0.65	0.56 - 0.69	0.58 - 0.63	0.52 - 0.63	0.53 - 0.63	0.54 - 0.59	0.55 - 0.74	0.5 - 0.6	0.58 - 0.59	0.51 - 0.69	0.61 - 0.7	
	Sigma: 5	0.65 - 0.71	0.51 - 0.54	0.6 - 0.65	0.62 - 0.72	0.52 - 0.66	0.55 - 0.61	0.54 - 0.66	0.54 - 0.62	0.5 - 0.61	0.59 - 0.74	0.51 - 0.82	0.63 - 0.7		
	Shape		0.61 - 0.7	0.53 - 0.59	0.55 - 0.59	0.62 - 0.66	0.58 - 0.6	0.56 - 0.58	0.52 - 0.6	0.55 - 0.68	0.53 - 0.67	0.57 - 0.62	0.52 - 0.69	0.52 - 0.54	
	WAV	HHH	0.58 - 0.66	0.52 - 0.58	0.61 - 0.68	0.64 - 0.69	0.55 - 0.64	0.55 - 0.56	0.51 - 0.56	0.59 - 0.71	0.5 - 0.57	0.57 - 0.62	0.58 - 0.64	0.51 - 0.57	
		HHL	0.58 - 0.66	0.51 - 0.57	0.59 - 0.6	0.64 - 0.7	0.55 - 0.64	0.55 - 0.58	0.57 - 0.64	0.54 - 0.65	0.5 - 0.67	0.61 - 0.65	0.51 - 0.7	0.59 - 0.66	
		HLH	0.55 - 0.66	0.51 - 0.55	0.56 - 0.64	0.59 - 0.68	0.55 - 0.57	0.54 - 0.59	0.58 - 0.65	0.57 - 0.65	0.5 - 0.56	0.57 - 0.61	0.55 - 0.67	0.53 - 0.61	
		HLL	0.6 - 0.7	0.5 - 0.61	0.57 - 0.6	0.6 - 0.63	0.55 - 0.57	0.58 - 0.63	0.52 - 0.63	0.54 - 0.6	0.5 - 0.66	0.6 - 0.63	0.5 - 0.69	0.54 - 0.65	
		LHH	0.65 - 0.66	0.54 - 0.61	0.63 - 0.66	0.58 - 0.69	0.54 - 0.58	0.55 - 0.59	0.53 - 0.62	0.55 - 0.62	0.5 - 0.56	0.58 - 0.62	0.53 - 0.7	0.54 - 0.62	
		LHL	0.64 - 0.68	0.51 - 0.56	0.58 - 0.66	0.59 - 0.7	0.54 - 0.58	0.54 - 0.57	0.51 - 0.66	0.54 - 0.72	0.5 - 0.68	0.59 - 0.65	0.53 - 0.7	0.52 - 0.66	
		LLH	0.59 - 0.66	0.51 - 0.61	0.52 - 0.63	0.63 - 0.7	0.54 - 0.57	0.57 - 0.62	0.53 - 0.58	0.54 - 0.66	0.5 - 0.68	0.57 - 0.61	0.52 - 0.64	0.58 - 0.61	
		LLL	0.62 - 0.66	0.51 - 0.58	0.59 - 0.64	0.56 - 0.6	0.54 - 0.61	0.51 - 0.58	0.55 - 0.65	0.56 - 0.7	0.5 - 0.61	0.55 - 0.59	0.52 - 0.68	0.58 - 0.61	
	CTD	Bin 64		0.58 - 0.64	0.51 - 0.53	0.56 - 0.64	0.62 - 0.72	0.56 - 0.7	0.54 - 0.55	0.53 - 0.6	0.54 - 0.72	0.5 - 0.57	0.56 - 0.61	0.53 - 0.72	0.6 - 0.63
		LOG	Sigma: 0.5	0.56 - 0.66	0.51 - 0.62	0.64 - 0.65	0.61 - 0.74	0.53 - 0.7	0.57 - 0.63	0.52 - 0.56	0.58 - 0.71	0.5 - 0.57	0.56 - 0.64	0.56 - 0.66	0.54 - 0.58
			Sigma: 1	0.59 - 0.64	0.5 - 0.54	0.62 - 0.63	0.53 - 0.74	0.55 - 0.6	0.51 - 0.57	0.53 - 0.62	0.55 - 0.62	0.5 - 0.61	0.56 - 0.64	0.54 - 0.64	0.54 - 0.57
Sigma: 1.5			0.58 - 0.64	0.5 - 0.53	0.66 - 0.68	0.61 - 0.74	0.55 - 0.63	0.52 - 0.56	0.51 - 0.59	0.56 - 0.62	0.5 - 0.58	0.56 - 0.64	0.53 - 0.68	0.57 - 0.58	
sigma: 2			0.63 - 0.73	0.5 - 0.55	0.61 - 0.64	0.67 - 0.74	0.55 - 0.62	0.57 - 0.58	0.51 - 0.55	0.6 - 0.68	0.5 - 0.62	0.58 - 0.64	0.55 - 0.66	0.53 - 0.57	
sigma: 2.5			0.62 - 0.66	0.5 - 0.55	0.56 - 0.61	0.6 - 0.74	0.55 - 0.62	0.56 - 0.58	0.51 - 0.66	0.55 - 0.76	0.5 - 0.62	0.56 - 0.64	0.54 - 0.65	0.53 - 0.59	
sigma: 3			0.59 - 0.64	0.51 - 0.55	0.55 - 0.58	0.62 - 0.74	0.55 - 0.61	0.55 - 0.6	0.51 - 0.66	0.56 - 0.62	0.5 - 0.62	0.57 - 0.64	0.56 - 0.71	0.61 - 0.65	
sigma: 3.5			0.62 - 0.68	0.5 - 0.56	0.57 - 0.61	0.62 - 0.74	0.55 - 0.59	0.53 - 0.55	0.56 - 0.74	0.57 - 0.69	0.5 - 0.61	0.58 - 0.64	0.55 - 0.7	0.56 - 0.63	
sigma: 4			0.56 - 0.67	0.5 - 0.52	0.66 - 0.69	0.61 - 0.74	0.55 - 0.59	0.51 - 0.55	0.56 - 0.66	0.58 - 0.72	0.5 - 0.65	0.57 - 0.66	0.56 - 0.81	0.62 - 0.63	
sigma: 4.5			0.61 - 0.71	0.51 - 0.56	0.57 - 0.59	0.67 - 0.74	0.55 - 0.64	0.55 - 0.6	0.58 - 0.65	0.56 - 0.65	0.5 - 0.69	0.57 - 0.64	0.51 - 0.76	0.54 - 0.69	
Sigma: 5		0.61 - 0.67	0.51 - 0.57	0.59 - 0.66	0.58 - 0.74	0.55 - 0.76	0.56 - 0.58	0.58 - 0.61	0.52 - 0.63	0.5 - 0.66	0.55 - 0.68	0.54 - 0.77	0.54 - 0.63		
Shape			0.59 - 0.62	0.51 - 0.55	0.54 - 0.55	0.6 - 0.7	0.6 - 0.61	0.54 - 0.54	0.56 - 0.56	0.57 - 0.64	0.53 - 0.66	0.56 - 0.67	0.54 - 0.69	0.56 - 0.56	
WAV		HHH	0.61 - 0.71	0.51 - 0.59	0.61 - 0.65	0.58 - 0.74	0.57 - 0.68	0.57 - 0.58	0.53 - 0.64	0.61 - 0.67	0.5 - 0.55	0.56 - 0.64	0.56 - 0.65	0.52 - 0.61	
		HHL	0.58 - 0.64	0.51 - 0.56	0.61 - 0.64	0.61 - 0.74	0.56 - 0.68	0.53 - 0.55	0.5 - 0.55	0.6 - 0.67	0.5 - 0.58	0.56 - 0.64	0.55 - 0.61	0.57 - 0.59	
		HLH	0.57 - 0.64	0.51 - 0.62	0.55 - 0.58	0.61 - 0.74	0.54 - 0.65	0.53 - 0.56	0.54 - 0.55	0.59 - 0.7	0.5 - 0.58	0.56 - 0.64	0.55 - 0.69	0.54 - 0.59	
		HLL	0.57 - 0.64	0.51 - 0.55	0.6 - 0.62	0.58 - 0.74	0.51 - 0.56	0.54 - 0.62	0.52 - 0.57	0.59 - 0.61	0.5 - 0.58	0.56 - 0.64	0.54 - 0.67	0.53 - 0.58	
		LHH	0.58 - 0.64	0.51 - 0.57	0.65 - 0.65	0.57 - 0.74	0.52 - 0.71	0.54 - 0.61	0.54 - 0.58	0.58 - 0.67	0.5 - 0.63	0.57 - 0.64	0.53 - 0.68	0.55 - 0.58	
		LHL	0.57 - 0.64	0.52 - 0.56	0.6 - 0.66	0.57 - 0.74	0.51 - 0.73	0.51 - 0.59	0.51 - 0.58	0.59 - 0.65	0.5 - 0.6	0.54 - 0.64	0.51 - 0.7	0.53 - 0.59	
		LLH	0.6 - 0.64	0.51 - 0.57	0.56 - 0.59	0.6 - 0.74	0.51 - 0.61	0.53 - 0.57	0.5 - 0.57	0.6 - 0.69	0.5 - 0.62	0.55 - 0.64	0.53 - 0.65	0.57 - 0.59	
		LLL	0.58 - 0.66	0.5 - 0.55	0.65 - 0.73	0.61 - 0.74	0.56 - 0.65	0.53 - 0.55	0.53 - 0.62	0.58 - 0.65	0.5 - 0.58	0.56 - 0.64	0.56 - 0.69	0.58 - 0.6	
PET		Bin 64		0.55 - 0.69	0.52 - 0.57	0.56 - 0.57	0.62 - 0.78	0.53 - 0.64	0.54 - 0.65	0.51 - 0.64	0.5 - 0.64	0.5 - 0.57	0.62 - 0.69	0.5 - 0.75	0.58 - 0.64
		LOG	Sigma: 0.5	0.58 - 0.68	0.5 - 0.52	0.62 - 0.68	0.62 - 0.74	0.54 - 0.68	0.54 - 0.6	0.54 - 0.62	0.51 - 0.65	0.5 - 0.67	0.56 - 0.61	0.5 - 0.81	0.61 - 0.67
			Sigma: 1	0.6 - 0.68	0.51 - 0.56	0.61 - 0.71	0.55 - 0.74	0.54 - 0.68	0.54 - 0.57	0.5 - 0.61	0.52 - 0.69	0.5 - 0.61	0.56 - 0.6	0.5 - 0.76	0.57 - 0.6
	Sigma: 1.5		0.57 - 0.68	0.5 - 0.54	0.64 - 0.7	0.61 - 0.74	0.51 - 0.65	0.54 - 0.56	0.5 - 0.64	0.5 - 0.64	0.5 - 0.6	0.56 - 0.61	0.5 - 0.68	0.55 - 0.6	
	sigma: 2		0.57 - 0.69	0.5 - 0.53	0.59 - 0.64	0.54 - 0.74	0.54 - 0.72	0.56 - 0.61	0.51 - 0.62	0.5 - 0.6	0.5 - 0.59	0.54 - 0.61	0.5 - 0.74	0.57 - 0.58	
	sigma: 2.5		0.64 - 0.69	0.51 - 0.56	0.57 - 0.68	0.6 - 0.74	0.54 - 0.61	0.56 - 0.61	0.51 - 0.61	0.5 - 0.61	0.5 - 0.58	0.6 - 0.65	0.5 - 0.67	0.53 - 0.58	
	sigma: 3		0.58 - 0.72	0.51 - 0.56	0.55 - 0.6	0.59 - 0.74	0.54 - 0.67	0.54 - 0.56	0.51 - 0.62	0.51 - 0.62	0.5 - 0.6	0.58 - 0.6	0.5 - 0.72	0.53 - 0.57	
	sigma: 3.5		0.62 - 0.71	0.51 - 0.56	0.54 - 0.56	0.58 - 0.74	0.54 - 0.62	0.55 - 0.61	0.51 - 0.64	0.5 - 0.65	0.5 - 0.57	0.6 - 0.69	0.5 - 0.75	0.56 - 0.62	
	sigma: 4		0.59 - 0.68	0.51 - 0.56	0.55 - 0.61	0.61 - 0.74	0.54 - 0.68	0.55 - 0.64	0.51 - 0.65	0.51 - 0.6	0.5 - 0.61	0.6 - 0.71	0.5 - 0.71	0.51 - 0.61	
	sigma: 4.5		0.57 - 0.68	0.52 - 0.56	0.6 - 0.64	0.56 - 0.74	0.55 - 0.64	0.54 - 0.67	0.51 - 0.78	0.51 - 0.75	0.5 - 0.56	0.6 - 0.69	0.5 - 0.79	0.53 - 0.62	
	Sigma: 5	0.67 - 0.69	0.51 - 0.56	0.55 - 0.6	0.58 - 0.74	0.55 - 0.65	0.54 - 0.6	0.51 - 0.68	0.5 - 0.64	0.5 - 0.56	0.62 - 0.67	0.5 - 0.75	0.55 - 0.67		
	Shape		0.59 - 0.72	0.51 - 0.61	0.57 - 0.62	0.6 - 0.66	0.57 - 0.57	0.56 - 0.57	0.53 - 0.53	0.54 - 0.64	0.51 - 0.61	0.6 - 0.67	0.53 - 0.66	0.53 - 0.53	
	WAV	HHH	0.65 - 0.71	0.51 - 0.61	0.6 - 0.81	0.64 - 0.74	0.5 - 0.62	0.52 - 0.61	0.5 - 0.61	0.51 - 0.62	0.5 - 0.58	0.56 - 0.61	0.5 - 0.76	0.57 - 0.62	
		HHL	0.62 - 0.74	0.51 - 0.56	0.62 - 0.7	0.65 - 0.74	0.51 - 0.69	0.52 - 0.58	0.51 - 0.6	0.5 - 0.59	0.5 - 0.57	0.59 - 0.62	0.5 - 0.63	0.53 - 0.62	
		HLH	0.61 - 0.68	0.51 - 0.56	0.57 - 0.64	0.66 - 0.74	0.52 - 0.64	0.57 - 0.61	0.51 - 0.65	0.51 - 0.67	0.5 - 0.61	0.61 - 0.67	0.5 - 0.67	0.53 - 0.65	
		HLL	0.61 - 0.68	0.51 - 0.56	0.59 - 0.63	0.61 - 0.74	0.54 - 0.61	0.54 - 0.58	0.52 - 0.57	0.5 - 0.61	0.5 - 0.65	0.55 - 0.61	0.5 - 0.69	0.56 - 0.61	
		LHH	0.6 - 0.68	0.52 - 0.56	0.61 - 0.74	0.63 - 0.74	0.51 - 0.55	0.54 - 0.54	0.51 - 0.56	0.5 - 0.65	0.5 - 0.6	0.56 - 0.61	0.5 - 0.68	0.53 - 0.56	
		LHL	0.61 - 0.68	0.51 - 0.53	0.58 - 0.6	0.58 - 0.74	0.52 - 0.6	0.56 - 0.6	0.51 - 0.63	0.5 - 0.56	0.5 - 0.58	0.57 - 0.61	0.5 - 0.61	0.57 - 0.6	
		LLH	0.61 - 0.68	0.51 - 0.56	0.58 - 0.62	0.58 - 0.74	0.52 - 0.72	0.53 - 0.67	0.51 - 0.6	0.51 - 0.56	0.5 - 0.6				

Supplemental Table 11. KRAS Classifier Mean+Sd different modality, preprocessing and setting

Modality	Preprocessed	Setting	AB	BAG	BNB	DT	GNB	KNN	LREG	MLP	QDA	RF	SGD	SVM	
CT	Bin 64		0.61 ± 0.02	0.55 ± 0.038	0.57 ± 0.019	0.63 ± 0.022	0.58 ± 0.04	0.54 ± 0.03	0.56 ± 0.017	0.6 ± 0.033	0.53 ± 0.048	0.57 ± 0.024	0.61 ± 0.08	0.59 ± 0.023	
	LOG	Sigma: 0.5	0.62 ± 0.044	0.51 ± 0.015	0.58 ± 0.012	0.61 ± 0.035	0.6 ± 0.049	0.55 ± 0.0076	0.58 ± 0.013	0.59 ± 0.018	0.54 ± 0.045	0.59 ± 0.014	0.62 ± 0.061	0.57 ± 0.015	
		Sigma: 1	0.66 ± 0.011	0.53 ± 0.019	0.6 ± 0.024	0.61 ± 0.023	0.56 ± 0.028	0.56 ± 0.0034	0.56 ± 0.028	0.59 ± 0.033	0.53 ± 0.039	0.57 ± 0.015	0.62 ± 0.066	0.54 ± 0.016	
		Sigma: 1.5	0.65 ± 0.0071	0.56 ± 0.027	0.59 ± 0.048	0.64 ± 0.038	0.56 ± 0.024	0.57 ± 0.019	0.55 ± 0.021	0.57 ± 0.023	0.56 ± 0.06	0.59 ± 0.023	0.61 ± 0.07	0.56 ± 0.017	
		sigma: 2	0.63 ± 0.032	0.54 ± 0.014	0.55 ± 0.032	0.61 ± 0	0.57 ± 0.018	0.57 ± 0.022	0.54 ± 0.026	0.58 ± 0.034	0.53 ± 0.039	0.59 ± 0.024	0.61 ± 0.077	0.56 ± 0.02	
		sigma: 2.5	0.64 ± 0.027	0.6 ± 0.011	0.58 ± 0.051	0.62 ± 0.026	0.57 ± 0.031	0.59 ± 0.028	0.56 ± 0.04	0.58 ± 0.018	0.54 ± 0.048	0.6 ± 0.018	0.59 ± 0.062	0.58 ± 0.018	
		sigma: 3	0.64 ± 0.019	0.55 ± 0.025	0.62 ± 0.019	0.6 ± 0.018	0.57 ± 0.027	0.54 ± 0.013	0.57 ± 0.027	0.64 ± 0.04	0.57 ± 0.072	0.59 ± 0.019	0.61 ± 0.08	0.66 ± 0.032	
		sigma: 3.5	0.64 ± 0.02	0.55 ± 0.025	0.64 ± 0.027	0.62 ± 0.025	0.57 ± 0.033	0.56 ± 0.01	0.68 ± 0.12	0.67 ± 0.062	0.54 ± 0.062	0.59 ± 0.035	0.66 ± 0.12	0.79 ± 0.021	
		sigma: 4	0.68 ± 0.033	0.54 ± 0.015	0.64 ± 0.044	0.65 ± 0.035	0.58 ± 0.04	0.64 ± 0.096	0.57 ± 0.056	0.64 ± 0.073	0.57 ± 0.084	0.59 ± 0.019	0.63 ± 0.094	0.63 ± 0.064	
		sigma: 4.5	0.64 ± 0.033	0.58 ± 0.049	0.61 ± 0.061	0.61 ± 0.021	0.6 ± 0.055	0.59 ± 0.042	0.55 ± 0.023	0.64 ± 0.07	0.53 ± 0.04	0.59 ± 0.0053	0.61 ± 0.066	0.66 ± 0.037	
		Sigma: 5	0.68 ± 0.029	0.52 ± 0.0095	0.61 ± 0.026	0.68 ± 0.035	0.58 ± 0.05	0.57 ± 0.023	0.57 ± 0.045	0.57 ± 0.036	0.54 ± 0.045	0.65 ± 0.053	0.64 ± 0.13	0.66 ± 0.034	
	Shape		0.64 ± 0.033	0.55 ± 0.022	0.57 ± 0.022	0.65 ± 0.016	0.59 ± 0.007	0.57 ± 0.0068	0.54 ± 0.029	0.62 ± 0.059	0.57 ± 0.052	0.61 ± 0.024	0.6 ± 0.063	0.53 ± 0.0072	
	WAV	HHH	0.64 ± 0.03	0.54 ± 0.023	0.65 ± 0.03	0.66 ± 0.023	0.57 ± 0.032	0.56 ± 0.0014	0.55 ± 0.019	0.65 ± 0.052	0.53 ± 0.03	0.59 ± 0.02	0.61 ± 0.027	0.53 ± 0.022	
		HHL	0.63 ± 0.034	0.54 ± 0.022	0.6 ± 0.0049	0.66 ± 0.023	0.6 ± 0.046	0.56 ± 0.012	0.59 ± 0.024	0.59 ± 0.039	0.57 ± 0.086	0.62 ± 0.016	0.61 ± 0.071	0.64 ± 0.025	
		HLH	0.63 ± 0.04	0.54 ± 0.018	0.59 ± 0.034	0.62 ± 0.035	0.55 ± 0.0092	0.58 ± 0.022	0.6 ± 0.025	0.61 ± 0.024	0.52 ± 0.024	0.59 ± 0.014	0.6 ± 0.043	0.58 ± 0.03	
		HLL	0.67 ± 0.039	0.56 ± 0.04	0.58 ± 0.013	0.62 ± 0.01	0.56 ± 0.0081	0.59 ± 0.019	0.59 ± 0.043	0.57 ± 0.021	0.56 ± 0.071	0.61 ± 0.013	0.6 ± 0.081	0.62 ± 0.046	
		LHH	0.65 ± 0.0059	0.56 ± 0.025	0.65 ± 0.011	0.63 ± 0.04	0.55 ± 0.017	0.57 ± 0.02	0.58 ± 0.044	0.6 ± 0.032	0.53 ± 0.03	0.6 ± 0.016	0.63 ± 0.059	0.6 ± 0.035	
		LHL	0.65 ± 0.014	0.54 ± 0.02	0.61 ± 0.036	0.65 ± 0.042	0.55 ± 0.016	0.56 ± 0.014	0.59 ± 0.05	0.62 ± 0.06	0.53 ± 0.07	0.62 ± 0.018	0.63 ± 0.069	0.61 ± 0.055	
		LLH	0.64 ± 0.027	0.56 ± 0.033	0.55 ± 0.053	0.66 ± 0.031	0.56 ± 0.011	0.59 ± 0.022	0.55 ± 0.025	0.6 ± 0.046	0.54 ± 0.068	0.6 ± 0.012	0.6 ± 0.055	0.6 ± 0.011	
		LLL	0.65 ± 0.015	0.55 ± 0.03	0.61 ± 0.021	0.58 ± 0.014	0.58 ± 0.033	0.56 ± 0.029	0.6 ± 0.04	0.61 ± 0.058	0.53 ± 0.042	0.57 ± 0.015	0.61 ± 0.069	0.59 ± 0.011	
		CTD	Bin 64		0.61 ± 0.023	0.52 ± 0.01	0.59 ± 0.04	0.66 ± 0.036	0.61 ± 0.051	0.54 ± 0.005	0.57 ± 0.029	0.61 ± 0.066	0.53 ± 0.037	0.58 ± 0.018	0.62 ± 0.069
	LOG	Sigma: 0.5	0.61 ± 0.034	0.56 ± 0.039	0.64 ± 0.0068	0.68 ± 0.049	0.62 ± 0.075	0.61 ± 0.028	0.54 ± 0.014	0.64 ± 0.054	0.53 ± 0.031	0.58 ± 0.027	0.62 ± 0.038	0.55 ± 0.016	
		Sigma: 1	0.61 ± 0.015	0.52 ± 0.014	0.63 ± 0.0025	0.64 ± 0.088	0.58 ± 0.017	0.53 ± 0.021	0.57 ± 0.033	0.59 ± 0.027	0.53 ± 0.043	0.58 ± 0.027	0.59 ± 0.041	0.56 ± 0.016	
		Sigma: 1.5	0.61 ± 0.019	0.52 ± 0.012	0.67 ± 0.0068	0.68 ± 0.048	0.59 ± 0.035	0.55 ± 0.022	0.53 ± 0.027	0.59 ± 0.027	0.53 ± 0.033	0.58 ± 0.029	0.58 ± 0.05	0.57 ± 0.0076	
sigma: 2		0.66 ± 0.038	0.53 ± 0.019	0.62 ± 0.011	0.7 ± 0.027	0.59 ± 0.033	0.58 ± 0.0069	0.54 ± 0.014	0.64 ± 0.027	0.53 ± 0.047	0.59 ± 0.021	0.61 ± 0.038	0.55 ± 0.014		
sigma: 2.5		0.63 ± 0.012	0.53 ± 0.019	0.57 ± 0.027	0.66 ± 0.059	0.59 ± 0.03	0.57 ± 0.0076	0.55 ± 0.054	0.62 ± 0.074	0.55 ± 0.058	0.59 ± 0.026	0.6 ± 0.044	0.56 ± 0.027		
sigma: 3		0.61 ± 0.016	0.53 ± 0.019	0.56 ± 0.014	0.67 ± 0.052	0.58 ± 0.02	0.57 ± 0.02	0.57 ± 0.052	0.58 ± 0.024	0.54 ± 0.051	0.59 ± 0.026	0.65 ± 0.054	0.62 ± 0.013		
sigma: 3.5		0.63 ± 0.024	0.53 ± 0.025	0.6 ± 0.02	0.67 ± 0.054	0.56 ± 0.016	0.54 ± 0.0059	0.6 ± 0.073	0.63 ± 0.056	0.53 ± 0.042	0.6 ± 0.02	0.63 ± 0.056	0.59 ± 0.038		
sigma: 4		0.62 ± 0.038	0.51 ± 0.0063	0.68 ± 0.018	0.67 ± 0.056	0.58 ± 0.02	0.53 ± 0.013	0.6 ± 0.035	0.64 ± 0.044	0.56 ± 0.065	0.61 ± 0.034	0.71 ± 0.08	0.62 ± 0.0055		
sigma: 4.5		0.66 ± 0.04	0.53 ± 0.016	0.58 ± 0.0078	0.69 ± 0.034	0.59 ± 0.034	0.57 ± 0.027	0.61 ± 0.023	0.61 ± 0.037	0.55 ± 0.074	0.59 ± 0.025	0.64 ± 0.086	0.62 ± 0.05		
Sigma: 5		0.64 ± 0.02	0.54 ± 0.02	0.64 ± 0.032	0.66 ± 0.066	0.6 ± 0.079	0.57 ± 0.011	0.59 ± 0.0098	0.6 ± 0.039	0.55 ± 0.064	0.61 ± 0.045	0.66 ± 0.083	0.61 ± 0.039		
Shape		0.61 ± 0.0091	0.53 ± 0.017	0.55 ± 0.0055	0.68 ± 0.04	0.6 ± 0.0055	0.54 ± 0	0.56 ± 0	0.59 ± 0.028	0.57 ± 0.044	0.59 ± 0.042	0.61 ± 0.049	0.56 ± 0		
WAV	HHH	0.64 ± 0.042	0.56 ± 0.03	0.64 ± 0.02	0.66 ± 0.066	0.63 ± 0.054	0.58 ± 0.0055	0.57 ± 0.041	0.64 ± 0.025	0.52 ± 0.021	0.59 ± 0.033	0.62 ± 0.03	0.58 ± 0.032		
	HHL	0.61 ± 0.019	0.53 ± 0.021	0.61 ± 0.014	0.67 ± 0.056	0.63 ± 0.061	0.54 ± 0.0082	0.53 ± 0.022	0.63 ± 0.025	0.53 ± 0.031	0.59 ± 0.033	0.6 ± 0.024	0.58 ± 0.011		
	HLH	0.61 ± 0.021	0.57 ± 0.038	0.56 ± 0.016	0.67 ± 0.052	0.6 ± 0.051	0.54 ± 0.014	0.54 ± 0.0067	0.64 ± 0.041	0.53 ± 0.031	0.58 ± 0.028	0.6 ± 0.054	0.57 ± 0.019		
	HLL	0.61 ± 0.021	0.53 ± 0.015	0.61 ± 0.01	0.66 ± 0.065	0.55 ± 0.018	0.59 ± 0.037	0.55 ± 0.017	0.6 ± 0.0084	0.53 ± 0.034	0.59 ± 0.027	0.62 ± 0.05	0.56 ± 0.019		
	LHH	0.61 ± 0.019	0.54 ± 0.017	0.65 ± 0	0.66 ± 0.07	0.62 ± 0.082	0.59 ± 0.029	0.56 ± 0.014	0.64 ± 0.032	0.54 ± 0.049	0.59 ± 0.024	0.61 ± 0.059	0.57 ± 0.01		
	LHL	0.61 ± 0.022	0.55 ± 0.016	0.62 ± 0.03	0.65 ± 0.069	0.62 ± 0.096	0.55 ± 0.032	0.54 ± 0.029	0.61 ± 0.024	0.53 ± 0.041	0.58 ± 0.033	0.61 ± 0.071	0.55 ± 0.023		
	LLH	0.61 ± 0.013	0.53 ± 0.021	0.58 ± 0.015	0.67 ± 0.055	0.56 ± 0.042	0.54 ± 0.012	0.54 ± 0.033	0.64 ± 0.033	0.55 ± 0.055	0.58 ± 0.03	0.59 ± 0.04	0.58 ± 0.01		
	LLL	0.61 ± 0.03	0.53 ± 0.018	0.67 ± 0.039	0.69 ± 0.047	0.61 ± 0.031	0.55 ± 0.014	0.58 ± 0.041	0.61 ± 0.028	0.53 ± 0.032	0.59 ± 0.027	0.64 ± 0.046	0.59 ± 0.0066		
	PET	Bin 64		0.62 ± 0.057	0.55 ± 0.019	0.56 ± 0.0082	0.68 ± 0.058	0.58 ± 0.042	0.56 ± 0.045	0.59 ± 0.057	0.6 ± 0.068	0.53 ± 0.031	0.66 ± 0.024	0.63 ± 0.11	0.62 ± 0.026
LOG	Sigma: 0.5	0.62 ± 0.043	0.51 ± 0.0079	0.65 ± 0.027	0.66 ± 0.043	0.6 ± 0.052	0.55 ± 0.022	0.57 ± 0.031	0.6 ± 0.059	0.56 ± 0.068	0.6 ± 0.019	0.66 ± 0.13	0.62 ± 0.022		
	Sigma: 1	0.64 ± 0.029	0.53 ± 0.017	0.68 ± 0.044	0.64 ± 0.068	0.59 ± 0.054	0.55 ± 0.011	0.57 ± 0.044	0.59 ± 0.065	0.54 ± 0.048	0.58 ± 0.015	0.66 ± 0.12	0.59 ± 0.0095		
	Sigma: 1.5	0.63 ± 0.036	0.52 ± 0.015	0.67 ± 0.026	0.66 ± 0.046	0.59 ± 0.057	0.55 ± 0.0063	0.58 ± 0.054	0.58 ± 0.062	0.53 ± 0.036	0.59 ± 0.018	0.59 ± 0.075	0.59 ± 0.019		
	sigma: 2	0.64 ± 0.053	0.52 ± 0.01	0.61 ± 0.023	0.64 ± 0.069	0.6 ± 0.068	0.57 ± 0.029	0.56 ± 0.041	0.56 ± 0.041	0.53 ± 0.036	0.59 ± 0.025	0.6 ± 0.094	0.57 ± 0.004		
	sigma: 2.5	0.65 ± 0.022	0.52 ± 0.018	0.63 ± 0.043	0.65 ± 0.052	0.58 ± 0.038	0.57 ± 0.023	0.55 ± 0.034	0.57 ± 0.043	0.53 ± 0.035	0.63 ± 0.022	0.59 ± 0.078	0.55 ± 0.02		
	sigma: 3	0.64 ± 0.051	0.54 ± 0.019	0.57 ± 0.031	0.66 ± 0.056	0.6 ± 0.053	0.55 ± 0.0089	0.55 ± 0.042	0.57 ± 0.045	0.53 ± 0.037	0.59 ± 0.0057	0.62 ± 0.1	0.56 ± 0.015		
	sigma: 3.5	0.66 ± 0.032	0.54 ± 0.021	0.55 ± 0.013	0.65 ± 0.055	0.59 ± 0.036	0.56 ± 0.024	0.59 ± 0.06	0.58 ± 0.057	0.54 ± 0.03	0.64 ± 0.036	0.64 ± 0.11	0.58 ± 0.024		
	sigma: 4	0.63 ± 0.032	0.54 ± 0.019	0.57 ± 0.028	0.66 ± 0.044	0.6 ± 0.054	0.56 ± 0.037	0.55 ± 0.056	0.55 ± 0.03	0.53 ± 0.043	0.63 ± 0.043	0.61 ± 0.094	0.58 ± 0.039		
	sigma: 4.5	0.64 ± 0.038	0.54 ± 0.02	0.62 ± 0.017	0.65 ± 0.063	0.6 ± 0.04	0.56 ± 0.05	0.6 ± 0.099	0.59 ± 0.086	0.53 ± 0.024	0.64 ± 0.033	0.62 ± 0.11	0.6 ± 0.042		
	Sigma: 5	0.68 ± 0.012	0.54 ± 0.021	0.58 ± 0.024	0.65 ± 0.053	0.59 ± 0.04	0.55 ± 0.022	0.58 ± 0.068	0.6 ± 0.067	0.52 ± 0.02	0.65 ± 0.016	0.63 ± 0.11	0.62 ± 0.043		
Shape		0.64 ± 0.051	0.56 ± 0.033	0.61 ± 0.026	0.63 ± 0.019	0.57 ± 0	0.56 ± 0.0069	0.53 ± 0.0019	0.58 ± 0.042	0.56 ± 0.033	0.63 ± 0.029	0.61 ± 0.044	0.53 ± 0.0022		
WAV	HHH	0.68 ± 0.026	0.57 ± 0.047	0.75 ± 0.1	0.68 ± 0.034	0.57 ± 0.047	0.56 ± 0.038	0.56 ± 0.044	0.58 ± 0.048	0.52 ± 0.033	0.59 ± 0.018	0.65 ± 0.12	0.6 ± 0.019		
	HHL	0.67 ± 0.04	0.53 ± 0.02	0.67 ± 0.036	0.68 ± 0.034	0.62 ± 0.083	0.53 ± 0.024	0.56 ± 0.036	0.56 ± 0.041	0.53 ± 0.036	0.61 ± 0.013	0.58 ± 0.062	0.58 ± 0.039		
	HLH	0.64 ± 0.028	0.52 ± 0.017	0.63 ± 0.036	0.68 ± 0.031	0.58 ± 0.046	0.58 ± 0.017	0.6 ± 0.068	0.59 ± 0.068	0.55 ± 0.052	0.63 ± 0.025	0.59 ± 0.073	0.61 ± 0.06		
	HLL	0.64 ± 0.025	0.52 ±												

Discussion

Accurate detection of EGFR/KRAS mutation status in NSCLCs allows improved selection of patients for effective therapeutic strategies (2). Although bio-techniques provide intra- and inter-tumor EGFR/KRAS mutation status in patients, they are not entirely feasible and suffer from limitations (38). In the present study, we developed a comprehensive radiomics framework that calculates radiomic features from low-dose CT, diagnostic CT, and PET images to predict EGFR and KRAS mutation status using univariate and multivariate ML algorithms in NSCLCs patients. In our radiogenomics study, several radiomics models were found as predictive for EGFR and KRAS mutation status, which may be utilized for non-invasive cost-effective assessment.

In recent years, radiogenomics—the correlation of imaging features with genomic parameters—have garnered significant interest, as non-invasive assessment frameworks to tailor therapy based on imaging biomarkers and underlying biological pathways (7, 8, 39). In our work, we observed a wide range of predictive model performances in different radiomics features, feature selection and classification algorithms, and for different imaging modalities. For EGFR and KRAS mutation status prediction univariate analysis of radiomics features and cross combination of imaging/validation /feature selection/Classification as multivariate analysis methods resulted in a wide range of performance (AUC: 0.5~0.83). The best predictive power of conventional PET parameters was achieved by SUV_{peak} (AUC: 0.69, P-value = 0.0002) and MTV (AUC: 0.55, P-value = 0.0011) for EGFR and KRAS. Univariate analysis of extracted radiomics features improved prediction power up to AUC: 0.75 (q-value: 0.003, Short Run Emphasis feature of GLRLM from LOG preprocessed image of PET with sigma value 1.5) and AUC: 0.71 (q-value 0.00005, The Large Dependence Low Gray Level Emphasis from GLDM in LOG preprocessed image of CTD sigma value 5) for EGFR and KRAS, respectively. Furthermore, the machine learning algorithm improved the prediction power of gene status up to AUC: 0.82 for EGFR (LOG preprocessed of PET image set with sigma 3 with VT feature selector and SGD classifier) and AUC: 0.83 for KRAS (CT image set with sigma 3.5 with SM feature selector and SGD classifier). Using radiomics feature and machine learning algorithm outperformed conventional methods in prediction of EGFR and KRAS gene status in NSCLC patients.

The role of feature selection is more prominent in cross-combination of feature selection and image-sets for EGFR mutation status prediction, and the Random Forest (RF) and Support Vector Machine (SVM) method has shown to yield a higher performance in EGFR and KRAS prediction. For classification, the selection of an appropriate algorithm has a more prominent role. Zhao *et al.* (40) developed a conventional radiomics model and 3D deep learning system to predict EGFR-mutant pulmonary adenocarcinoma in CT images, the best AUC were reported by 0.758 using deep learning, however, in present study the of EGFR status achieved to AUC: 0.82 which shows the role of different combination of image set, feature selector and classifier. Parmer *et al.* (24) compared 14 feature selection and 12 classification methods in terms of their performance and stability, and found that the Wilcoxon rank-sum test and random forest (RF) had the highest prognostic performance for feature selection and classification, respectively, while demonstrating high stability against data perturbation. Moreover, Zhang *et al.* (41) built a 54-cross combination ML algorithms framework including six feature selections and nine classification methods for survival prediction of advanced nasopharyngeal carcinoma. They found that using RF for both feature selection and classification had the highest prognostic performance, followed by RF+AB and sure independence screening (SIS) + linear support vector machines (SVMs). Abdollahi *et al.* (42) studied different combination of feature selector and classifier for prediction of intensity-modulated radiation therapy response, Gleason score and stage in prostate cancer and they showed different task of classification can be achieved by different combination of feature selector and classifiers. In addition, we compared the predictive power of conventional PET parameters and radiomic features. In our radiogenomics study, we showed that radiomics features and machine learning based radiomic models are more predictive than conventional parameters.

Radiomics/genomics studies suffer from several challenges, and a robust framework for clinical decision making is highly desired (7, 43, 44). As an approved guideline, standardization efforts (IBSI in particular)(36) have sought to address the challenge of reproducing and validating reported findings by comparing and standardizing definitions and implementation of several image-feature sets between participating institutions. The IBSI also provides image biomarker nomenclature and definitions, benchmark data sets, and benchmark values to verify image processing and image biomarker calculations, as well as reporting guidelines, for high-throughput image analysis (36). In this study, we followed the IBSI protocol to resample all images to a constant voxel size to reduce feature sensitivity to variable image generation parameters (36), also ensuring the robustness and reproducibility of the study (45). Having a dataset of images generated using varying data acquisition (46, 47), reconstruction (48), processing (49) and segmentation (50) results produce inconsistencies in feature evaluation.

The limited size of a dataset is a limiting factor in radiomics studies, as the analysis with many radiomic features and only a few imaging data is prone to over-fitting. We tuned our models using 10-fold cross-validation to reduce the sensitivity of our results to input data and repeated for 20 times to make our results more reliable and all model evaluation were performed on 63 patients' independent validation set which never used in training process. Our presented results were developed using a standardized radiomics analysis workflow that showed strong and significant predictability of gene mutations. Specifically, in the feature selection phase, we used ML algorithms to reduce the number of features in order to decrease the dimensionality and reduce over-fitting.

Conclusion

We demonstrated that EGFR and KRAS mutation status in NSCLC patients can be predicted via a non-invasive and reliable radiomics analysis. We evaluated several different machine learning methods to find optimal methods for radiogenomics analyses. Overall, we showed that radiomic features extracted from different imaging modalities (alone or in combination) could be used for successful prediction of EGFR and KRAS mutation status, having more predictive power than conventional imaging biomarkers.

Reference:

1. Roberts PJ, Stinchcombe TE, Der CJ, Socinski MA. Personalized medicine in non–small-cell lung cancer: is KRAS a useful marker in selecting patients for epidermal growth factor receptor–targeted therapy? *Journal of clinical oncology*. 2010;28(31):4769-77.
2. Linardou H, Dahabreh IJ, Kanaloupiti D, Siannis F, Bafaloukos D, Kosmidis P, et al. Assessment of somatic k-RAS mutations as a mechanism associated with resistance to EGFR-targeted agents: a systematic review and meta-analysis of studies in advanced non-small-cell lung cancer and metastatic colorectal cancer. *The lancet oncology*. 2008;9(10):962-72.
3. Eberhard DA, Johnson BE, Amler LC, Goddard AD, Heldens SL, Herbst RS, et al. Mutations in the epidermal growth factor receptor and in KRAS are predictive and prognostic indicators in patients with non–small-cell lung cancer treated with chemotherapy alone and in combination with erlotinib. *J Clin Oncol*. 2005;23(25):5900-9.
4. Riely GJ, Marks J, Pao W. KRAS mutations in non–small cell lung cancer. *Proceedings of the American Thoracic Society*. 2009;6(2):201-5.
5. Mak R, Hermann G, Aerts H, Chen A, Baldini E, Kozono D, et al. Outcomes by EGFR, KRAS and ALK Genotype After Combined Modality Therapy for Locally Advanced Non-Small Cell Lung Cancer. *International Journal of Radiation Oncology• Biology• Physics*. 2016;96(2):S156.
6. Sun X, Xiao Z, Chen G, Han Z, Liu Y, Zhang C, et al. A PET imaging approach for determining EGFR mutation status for improved lung cancer patient management. *Sci Transl Med*. 2018;10(431):eaan8840.
7. Wu J, Tha KK, Xing L, Li R. Radiomics and radiogenomics for precision radiotherapy. *Journal of radiation research*. 2018;59(suppl_1):i25-i31.
8. Thawani R, McLane M, Beig N, Ghose S, Prasanna P, Velcheti V, et al. Radiomics and Radiogenomics in Lung Cancer: A Review for the Clinician. *Lung Cancer*. 2017.
9. Incoronato M, Aiello M, Infante T, Cavaliere C, Grimaldi AM, Mirabelli P, et al. Radiogenomic analysis of oncological data: a technical survey. *International journal of molecular sciences*. 2017;18(4):805.
10. Lambin P, Rios-Velazquez E, Leijenaar R, Carvalho S, van Stiphout RG, Granton P, et al. Radiomics: extracting more information from medical images using advanced feature analysis. *European journal of cancer*. 2012;48(4):441-6.
11. Gillies RJ, Kinahan PE, Hricak H. Radiomics: images are more than pictures, they are data. *Radiology*. 2015;278(2):563-77.
12. Hatt M, Tixier F, Pierce L, Kinahan PE, Le Rest CC, Visvikis D. Characterization of PET/CT images using texture analysis: the past, the present... any future? *European journal of nuclear medicine and molecular imaging*. 2017;44(1):151-65.
13. Beik J, Shiran MB, Abed Z, Shiri I, Ghadimi-Daresajini A, Farkhondeh F, et al. Gold nanoparticle-induced sonosensitization enhances the antitumor activity of ultrasound in colon tumor-bearing mice. *Medical physics*. 2018;45(9):4306-14.
14. Mirrahimi M, Abed Z, Beik J, Shiri I, Dezfuli AS, Mahabadi VP, et al. A thermo-responsive alginate nanogel platform co-loaded with gold nanoparticles and cisplatin for combined cancer chemo-photothermal therapy. *Pharmacological research*. 2019.
15. Zinn PO, Majadan B, Sathyan P, Singh SK, Majumder S, Jolesz FA, et al. Radiogenomic mapping of edema/cellular invasion MRI-phenotypes in glioblastoma multiforme. *PloS one*. 2011;6(10):e25451.
16. Kang J, Rancati T, Lee S, Oh JH, Kerns SL, Scott JG, et al. Machine learning and radiogenomics: lessons learned and future directions. *Frontiers in oncology*. 2018;8:228.
17. Abdollahi H, Shiri I, Heydari M. Medical Imaging Technologists in Radiomics Era: An Alice in Wonderland Problem. *Iranian journal of public health*. 2019;48(1):184.
18. Abdollahi H, Mahdavi SR, Shiri I, Mofid B, Bakhshandeh M, Rahmani K. Magnetic resonance imaging radiomic feature analysis of radiation-induced femoral head changes in prostate cancer radiotherapy. *Journal of cancer research and therapeutics*. 2019;15(8):11.
19. Liu Y, Kim J, Balagurunathan Y, Li Q, Garcia AL, Stringfield O, et al. Radiomic features are associated with EGFR mutation status in lung adenocarcinomas. *Clin Lung Cancer*. 2016;17(5):441-8. e6.
20. Lin YC, Lin G, Hong JH, Lin YP, Chen FH, Ng SH, et al. Diffusion radiomics analysis of intratumoral heterogeneity in a murine prostate cancer model following radiotherapy: Pixelwise correlation with histology. *J Magn Reson Imaging*. 2017;46(2):483-9.
21. Velazquez ER, Parmar C, Liu Y, Coroller TP, Cruz G, Stringfield O, et al. Somatic mutations drive distinct imaging phenotypes in lung cancer. *Cancer Res*. 2017;77(14):3922-30.
22. Zhang L, Chen B, Liu X, Song J, Fang M, Hu C, et al. Quantitative Biomarkers for Prediction of Epidermal Growth Factor Receptor Mutation in Non-Small Cell Lung Cancer. *Transl Oncol*. 2018;11(1):94-101.
23. Lambin P, Leijenaar RT, Deist TM, Peerlings J, de Jong EE, van Timmeren J, et al. Radiomics: the bridge between medical imaging and personalized medicine. *Nature Reviews Clinical Oncology*. 2017;14(12):749.
24. Parmar C, Grossmann P, Bussink J, Lambin P, Aerts HJ. Machine learning methods for quantitative radiomic biomarkers. *Scientific reports*. 2015;5:13087.
25. Parmar C, Grossmann P, Rietveld D, Rietbergen MM, Lambin P, Aerts HJ. Radiomic machine-learning classifiers for prognostic biomarkers of head and neck cancer. *Frontiers in oncology*. 2015;5:272.
26. Leger S, Zwanenburg A, Pilz K, Lohaus F, Linge A, Zöphel K, et al. A comparative study of machine learning methods for time-to-event survival data for radiomics risk modelling. *Scientific reports*. 2017;7(1):13206.
27. Deist TM, Dankers FJ, Valdes G, Wijsman R, Hsu IC, Oberije C, et al. Machine learning algorithms for outcome prediction in (chemo) radiotherapy: An empirical comparison of classifiers. *Medical physics*. 2018.
28. Abdollahi H, Mostafaei S, Cheraghi S, Shiri I, Mahdavi SR, Kazemnejad A. Cochlea CT radiomics predicts chemoradiotherapy induced sensorineural hearing loss in head and neck cancer patients: a machine learning and multi-variable modelling study. *Physica Medica*. 2018;45:192-7.
29. Prior FW, Clark K, Commean P, Freymann J, Jaffe C, Kirby J, et al., editors. TCIA: an information resource to enable open science. *Conf Proc IEEE Eng Med Biol Soc*; 2013: IEEE.
30. Gevaert O, Xu J, Hoang CD, Leung AN, Xu Y, Quon A, et al. Non–small cell lung cancer: identifying prognostic imaging biomarkers by leveraging public gene expression microarray data—methods and preliminary results. *Radiology*. 2012;264(2):387-96.
31. Clark K, Vendt B, Smith K, Freymann J, Kirby J, Koppel P, et al. The Cancer Imaging Archive (TCIA): maintaining and operating a public information repository. *Journal of digital imaging*. 2013;26(6):1045-57.
32. Bakr S, Gevaert O, Echegaray S, Ayers K, Zhou M, Shafiq M, et al. A radiogenomic dataset of non-small cell lung cancer. *Scientific Data*. 2018;5:180202.
33. Rosset A, Spadola L, Ratib O. OsiriX: an open-source software for navigating in multidimensional DICOM images. *Journal of digital imaging*. 2004;17(3):205-16.
34. Pieper S, Halle M, Kikinis R, editors. 3D Slicer. *Biomedical Imaging: Nano to Macro, 2004 IEEE International Symposium on*; 2004: IEEE.
35. van Griethuysen JJ, Fedorov A, Parmar C, Hosny A, Aucoin N, Narayan V, et al. Computational Radiomics System to Decode the Radiographic Phenotype. *Cancer research*. 2017;77(21):e104-e7.

36. Zwanenburg A, Leger S, Vallières M, Löck S. Image biomarker standardisation initiative. arXiv preprint arXiv:161207003. 2016.
37. Pedregosa F, Varoquaux G, Gramfort A, Michel V, Thirion B, Grisel O, et al. Scikit-learn: Machine learning in Python. *Journal of machine learning research*. 2011;12(Oct):2825-30.
38. Simon R, Roychowdhury S. Implementing personalized cancer genomics in clinical trials. *Nature reviews Drug discovery*. 2013;12(5):358.
39. Shiri I, Maleki H, Hajianfar G, Abdollahi H, Ashrafinia S, Hatt M, et al. PET/CT Radiomic Sequencer for Prediction of EGFR and KRAS Mutation Status in NSCLC Patients. arXiv e-prints [Internet]. 2019 June 01, 2019. Available from: <https://ui.adsabs.harvard.edu/abs/2019arXiv190606623S>.
40. Zhao W, Yang J, Ni B, Bi D, Sun Y, Xu M, et al. Toward automatic prediction of EGFR mutation status in pulmonary adenocarcinoma with 3D deep learning. *Cancer medicine*. 2019.
41. Zhang B, He X, Ouyang F, Gu D, Dong Y, Zhang L, et al. Radiomic machine-learning classifiers for prognostic biomarkers of advanced nasopharyngeal carcinoma. *Cancer letters*. 2017;403:21-7.
42. Abdollahi H, Mofid B, Shiri I, Razzaghdoost A, Saadipoor A, Mahdavi A, et al. Machine learning-based radiomic models to predict intensity-modulated radiation therapy response, Gleason score and stage in prostate cancer. *La radiologia medica*. 2019:1-13.
43. Reuzé S, Schernberg A, Orhac F, Sun R, Chargari C, Dercle L, et al. Radiomics in Nuclear Medicine applied to Radiation Therapy: methods, pitfalls and challenges. *International Journal of Radiation Oncology• Biology• Physics*. 2018.
44. Aerts HJ. The potential of radiomic-based phenotyping in precision medicine: a review. *JAMA oncology*. 2016;2(12):1636-42.
45. Pinker K, Shitano F, Sala E, Do RK, Young RJ, Wibmer AG, et al. Background, current role, and potential applications of radiogenomics. *Journal of Magnetic Resonance Imaging*. 2018;47(3):604-20.
46. Leijenaar RT, Carvalho S, Velazquez ER, Van Elmpt WJ, Parmar C, Hoekstra OS, et al. Stability of FDG-PET Radiomics features: an integrated analysis of test-retest and inter-observer variability. *Acta oncologica*. 2013;52(7):1391-7.
47. Shiri I, Abdollahi H, Shaysteh S, Mahdavi SR. Test-Retest Reproducibility and Robustness Analysis of Recurrent Glioblastoma MRI Radiomics Texture Features. *Iranian Journal of Radiology*. 2017(5).
48. Shiri I, Rahmim A, Ghaffarian P, Geramifar P, Abdollahi H, Bitarafan-Rajabi A. The impact of image reconstruction settings on 18F-FDG PET radiomic features: multi-scanner phantom and patient studies. *European radiology*. 2017;27(11):4498-509.
49. Shiri I, Ghafarian P, Geramifar P, Leung KH-Y, Ghelichoghli M, Oveisi M, et al. Direct attenuation correction of brain PET images using only emission data via a deep convolutional encoder-decoder (Deep-DAC). *European radiology*. 2019:1-13.
50. Lu L, Lv W, Jiang J, Ma J, Feng Q, Rahmim A, et al. Robustness of Radiomic Features in [11 C] Choline and [18 F] FDG PET/CT Imaging of Nasopharyngeal Carcinoma: Impact of Segmentation and Discretization. *Molecular Imaging and Biology*. 2016;18(6):935-45.