

Bending modes metrology in the 14-15 μm region

M. Lamperti¹, R. Gotti¹, D. Gatti¹, M. K. Shakfa², E. Cané³, F. Tamassia³, P. Schunemann⁴, P. Laporta¹, A. Farooq^{2*}, and M. Marangoni^{1*}

¹*Dipartimento di Fisica - Politecnico di Milano and IFN-CNR, Via Gaetano Previati 1/C, 23900 Lecco, Italy*

²*King Abdullah University for Science and Technology, Clean Combustion Research Center, Thuwal 23955, Saudi Arabia*

³*Università di Bologna, Dipartimento di Chimica Industriale, Viale del Risorgimento 4, 40136 Bologna, Italy*

⁴*BAE Systems, Inc., MER15-1813, P.O. Box 868, Nashua, New Hampshire 03061-0868, USA*

**aamir.farooq@kaust.edu.sa, marco.marangoni@polimi.it*

Frequency combs have triggered an impressive evolution of optical metrology across diverse regions of the electromagnetic spectrum, from the ultraviolet to the terahertz frequencies^{1,2}. An unexplored territory, however, remains in the region of vibrational bending modes, mostly due to the lack of single-mode lasers in the long-wavelength (LW) part of the mid-infrared (MIR) spectrum³. We fill this gap through a purely MIR-based nonlinear laser source with tunability from 12.1 to 14.8 μm , optical power up to 110 μW , MHz-level linewidth and comb calibration. This enables the first example of bending modes metrology in this region, with the assessment of several CO₂-based frequency benchmarks with uncertainties down to 30 kHz, and the accurate study of the ν_{11} band of benzene, which is a significant testbed for the resolution of the spectrometer. These achievements pave the way for LW-MIR metrology⁴, rotationally-resolved studies⁵ and astronomic observations⁶ of large molecules, such as aromatic hydrocarbons.

The lack of continuously tunable single-mode lasers in the region from 13 to 20 μm has been a stumbling block for optical metrology and high-resolution spectroscopy in the region of vibrational bending modes⁷. Cryogenically cooled lead-salt diode lasers have almost disappeared in recent years⁸. They have been replaced by room-temperature distributed-feedback quantum-cascade-lasers (DFB-QCLs) that offer portability, stable operation and mode-hop-free tuning⁹. However, they are commercially available only below 13 μm ¹⁰, whereas the emerging technology of InAs-based long-wavelength QCLs is still under development and not available for applications¹¹. An alternative approach is difference frequency generation (DFG) from continuous-wave (cw) Ti:sapphire^{12,13} or diode¹⁴ lasers, but this is accompanied by extremely low optical powers (10-100 nW range) that impair the acquisition of absorption spectra at high signal-to-noise ratio (SNR). The limitations of laser technology have established Fourier-Transform spectroscopy driven by incoherent light sources as a gold standard for the LW-MIR region, at the price of low resolution (0.0007 cm^{-1} in the best cases¹⁵) and of no absolute calibration for the frequency axis. A powerful solution to both issues is direct comb spectroscopy¹⁶, which has recently conquered the LW-MIR region by a dual-comb approach¹⁷ that enabled snapshots of entire bands from 6.7 to 16.7 μm at high temporal and spectral resolution¹⁸. This achievement does not remove, however, the need for cw probe lasers to observe single lines with high accuracy^{4,7,19} and extreme speed²⁰ in a region where the Doppler width and the typical line densities are consistently below the spacing between adjacent comb modes. Furthermore, to the best of our knowledge, optical metrology has not been demonstrated so far by direct comb spectroscopy in the LW-MIR region.

Our spectrometer is schematically shown in Fig. 1a. The laser source is based on the DFG process between two MIR sources, namely a cw DFB-QCL and a CO₂ laser, in an orientation-patterned gallium arsenide (OP-GaAs) crystal²¹. The DFG laser provides highly coherent radiation between 12.1 and 14.8 μm with optical power up to 110 μW and a 2.2 MHz linewidth. The large spectral range comes from the coarse wavelength adjustment of the CO₂ laser (from 9.23 to 10.86 μm) combined with the fine temperature-based tuning of the QCL (over ~4 cm⁻¹ around 5.69 μm) and the fan-out structure of the poling periods of the OP-GaAs crystal (see Methods for details). The optical power is 2-3 orders of magnitude better than previous DFG sources starting from cw Ti:sapphire lasers with GaSe crystals^{12,13}. The obtained power agrees within 15% of our calculations using 63 pm/V as the effective nonlinear coefficient. Such an excellent agreement reflects the quality of the crystal and a very favourable MIR-MIR interaction that ensures similar diffraction conditions for the interacting beams and thus optimal spatial overlap over a long crystal. The absolute frequency calibration of the DFG radiation comes from the referencing of both QCL and CO₂ lasers to a 100 MHz comb *via* sum frequency generation²² (see Methods for details). As sketched in Fig. 1b, the CO₂ laser remains offset-locked to the nearest comb mode during the measurements, whereas the frequency of the temperature-tuned DFB-QCL is tracked against the comb by real-time fast Fourier transform (FFT) analysis and barycentre calculation of its beat note. Synchronous acquisition of the beat note and gas transmission allows spectra to be straightforwardly calibrated.

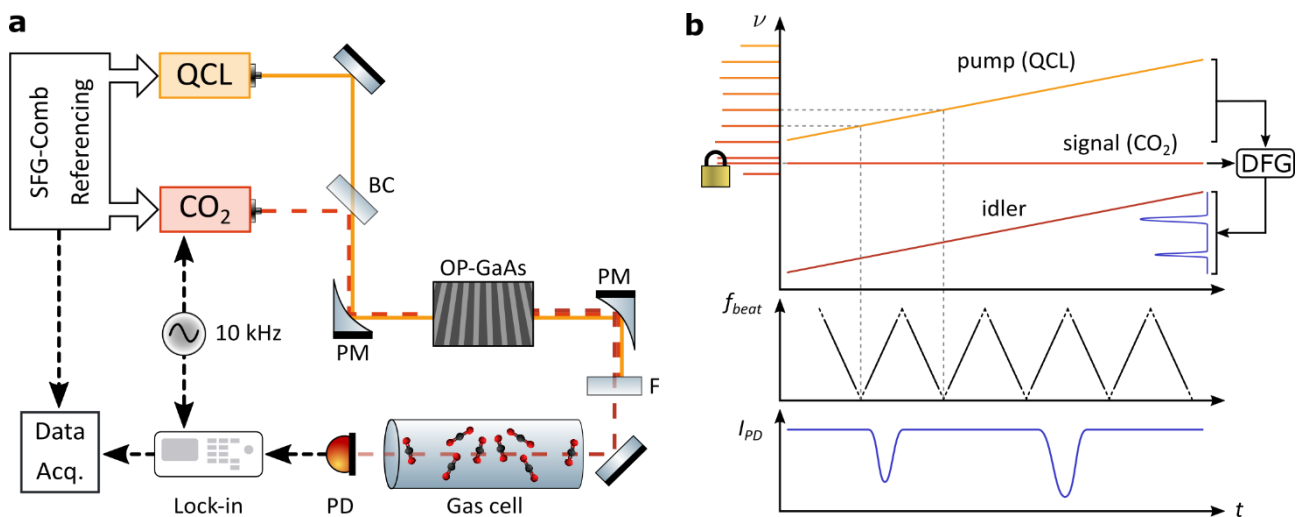


Fig. 1 | Comb-referenced nonlinear laser source for the LW-MIR. a. Experimental lay-out. BC: beam combiner; PM: parabolic mirror; F: low-pass filter; PD: liquid-nitrogen cooled photodiode; SFG: sum frequency generation. **b.** Pictorial view of the DFG process between the comb-locked CO₂ laser and the QCL. The synchronous acquisition of the QCL-comb beat note (f_{beat}) and of the idler absorption (I_{PD}) allow absolute calibration of the spectra over the measurement time. For graphical simplicity pump and signal are shown to beat directly with a MIR comb, but in practice this happens through the interposition of a sum frequency generation process (see Methods for details).

Figure 2a shows an example of comb-calibrated transmission spectrum for some intense lines of the ν_2 band of CO₂ in the 675-689 cm⁻¹ (14.5-14.8 μm) range. Lines are well isolated and in excellent qualitative agreement with the HITRAN²³ simulations.

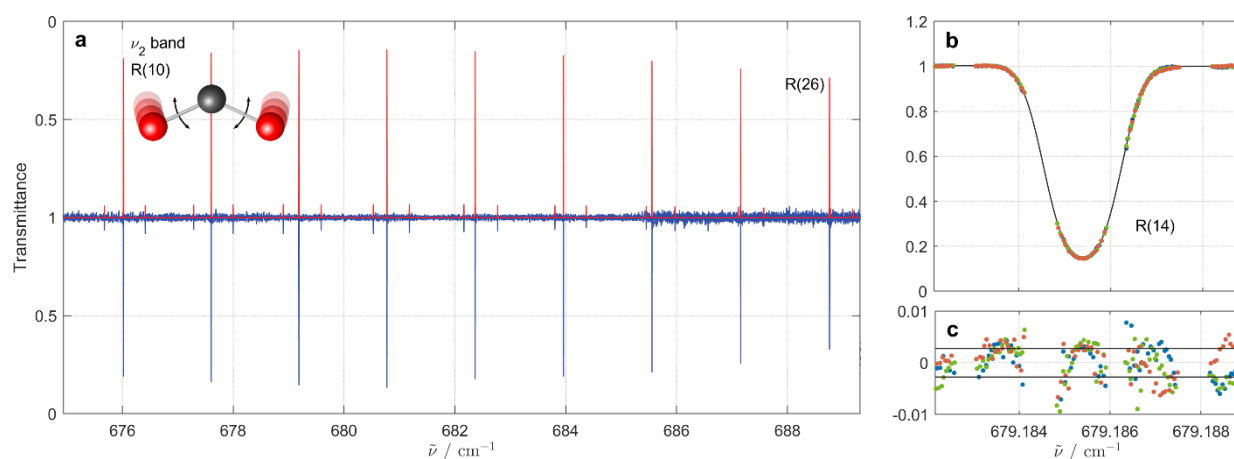


Fig. 2 | Carbon dioxide spectroscopy. **a.** Absorption spectrum of the ν_2 band of CO_2 (blue) as compared to HITRAN simulations (red). Experimental conditions: interaction length = 0.67 m, pressure = 0.93 Pa, room temperature. The scattering of baseline points is due to etalons that have been fitted out in the proximity of lines. **b.** Zoomed-in view of the R(14) line, with independent measurements of the same line (coloured dots), and fit with a Voigt profile (black line). Each measurement is composed of 170 points acquired over 2 seconds. **c.** Residuals from the Voigt fitting of each individual measurement, showing the presence of a periodic structure due to parasitic etalons in the setup. The black lines delimit the region of ± 1 standard deviation of the residuals (0.3 % of the line amplitude).

Figure 2b shows three independent experimental traces of the R(14) line, reported in different colors to highlight their reproducibility over the two measurement axes. The presence of spectral gaps is due to the missed calibration that occurs when the comb-QCL beat note approaches either DC or half the spacing between neighboring comb modes (50 MHz, see Methods for details). With an additional spectral acquisition at a slightly different comb repetition rate, the gaps would disappear. They do not, however, represent an issue when fitting the data (black line in Fig. 2b): the residuals in Fig. 2c refer to a Voigt model and give an SNR of approximately 300 with a measurement time of 2 seconds. This translates to a relatively small statistical uncertainty of about 140 kHz on the inferred line center frequency, mostly due to the reduced Doppler broadening (30 MHz level) given by the LW-MIR range and to the negligible collisional broadening at the low pressure used here.

For each measured line, Fig. 3 shows the difference between its absolute line-centre frequency and the value in the HITRAN database, together with an error bar dominated by statistical uncertainties (see Methods for details). Some lines have been intentionally measured several times to improve their line-centre uncertainty and to provide frequency standards in a region where absolute calibration has not been reported so far, to the best of our knowledge. The weighted deviation from HITRAN, averaged over 24 transitions, amounts to a remarkably low 0.34 MHz, which is consistent with the declared HITRAN uncertainty (0.3-3 MHz). The line centres reported for each line in the Supplementary Table 1 offer a metrological benchmark to calibrate past and future spectra in the $670\text{-}720\text{ cm}^{-1}$ (14-15 μm) region.

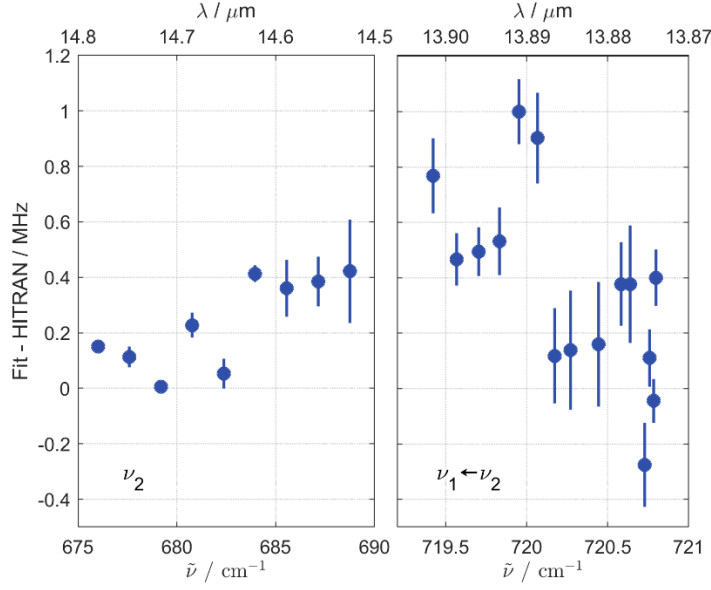


Fig. 3 | Carbon dioxide: experimental line centres vs HITRAN values. Difference between experimental and HITRAN line centres for several lines of the R branch of the ν_2 band (left panel) and of the Q branch of the $\nu_1 \leftarrow \nu_2$ band (right panel). The variation of error bars comes from different number of measurements per line.

In Fig. 4, we report the experimental spectrum of the R branch of the ν_{11} band of benzene over 675-688 cm^{-1} . This intense, parallel band is an interesting testbed for the spectral resolution and frequency repeatability of the spectrometer because of the large density of rotational lines. Apart from transitions with quantum number $K = 0, 1$ and 2, whose separation is smaller than the Doppler width, all K lines are well resolved. This is shown in the insets in Fig. 4 and is particularly evident at large J . Lines from three hot bands, namely $\nu_{11} + \nu_{16} - \nu_{16}$, $\nu_{11} + \nu_6 - \nu_6$ and $2\nu_{11} - \nu_{11}$, are also present in the recorded spectral region, but the analysis was restricted to the 729 ro-vibration transitions of the ν_{11} band, with J up to 38 and K up to 36, since these are mostly isolated, medium intensity lines. We adopted a weighted least-squares fitting procedure to determine highly accurate term values of the excited ro-vibration levels, according to the equation:

$$E_v(J, k) = E_v^0 + B_v[J(J+1) - k^2] + C_v k^2 - D_{v,J}[J(J+1)]^2 - D_{v,JK}J(J+1)k^2 - D_{v,K}k^4 + H_{v,J}[J(J+1)]^3 + H_{v,JK}[J(J+1)]^2 k^2 + H_{v,KJ}[J(J+1)]k^4 + H_{v,K}k^6 \quad (1)$$

where v stands for $\nu_{11} = 1$ and $k = \pm K$. The ground state term values were derived from Ref. 24 and Eq. (1) with $v = 0$, $E_v^0 = 0$ and keeping all H centrifugal distortion parameters fixed to zero. A unitary weight was attributed to the isolated lines, whose wavenumber precision is estimated to be $1 \times 10^{-5} \text{ cm}^{-1}$. The weights of blended or weak lines were reduced to 0.01, and, in addition, for lines with multiple assignments, the weight was further reduced by a factor equal to the number of overlapped transitions. The parameters obtained from the best fit are listed in Table 1, along with those of the ground state for comparison.

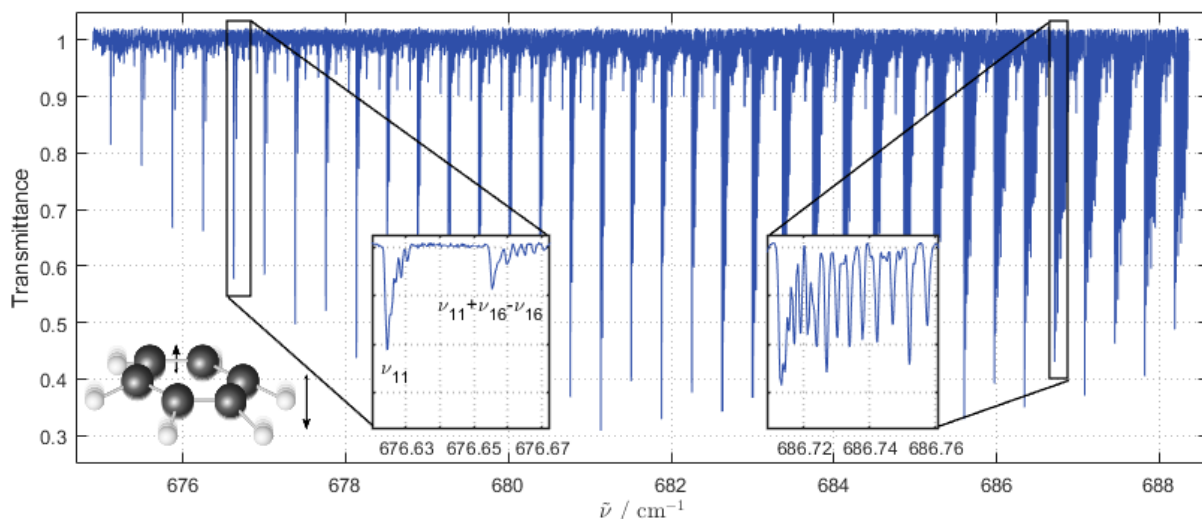


Fig. 4 | Absorption spectrum of the R-branch of the ν_{11} band of benzene. Experimentally measured transmission at the output of a 0.67 m long cell at a pressure of 13.3 Pa and room temperature. The spectral features have a typical shape, with the K -structure degrading regularly towards higher wavenumbers and an apparent intensity alternation of subsequent lines. The insets show zoomed-in views of $R_K(6)$ and $R_K(33)$ features (left and right inset, respectively): for the latter, the K -splitting is much more manifest together with intensity alternation of K components, owing to spin statistics, that varies as 10:11:9:14, for $K = 6p, 6p\pm 1, 6p\pm 2, 6p\pm 3$, respectively, ($p = 0, 1, 2, 3, \dots$)²⁵.

All parameters reported in Table 1 are statistically well determined. The uncertainty of the band origin (E_v^0) has been reduced by about one third with respect to its value in Ref. 26. The B, C, D_J, D_{JK} , and D_K constants are consistent and close to those of the ground state²⁴. The values of B, D_J , and D_{JK} are more accurate than those in Ref. 26 by one order of magnitude. The values of C, D_K , and H_J constants of $\nu_{11} = 1$ are determined here for the first time. The standard deviation of the fit, $5.1 \times 10^{-5} \text{ cm}^{-1}$ (1.5 MHz), is one order of magnitude smaller than in the past literature for the same band²⁶. At the same time, it is about 10 times higher than for single CO_2 lines due to a combination of factors, namely, the many partially overlapped lines and the difficulty to access intensity baseline due to the congested spectrum. The absolute calibration of the frequency axis allows, for the first time, the stacking of repeated measurements and the determination of absolute line positions at these wavelengths.

Table 1. Spectroscopic parameters (cm^{-1}) for the ground state and for the $\nu_{11} = 1$ excited state of benzene^a

		Ground state	$\nu_{11} = 1$
E_v^0			673.9751463(121)
B		0.1897727	0.1896373408(614)
C		0.0948863	0.0949217259(632)
D_J	$\times 10^8$	3.34	3.41413(847)
D_{JK}	$\times 10^8$	-6.57	-6.54650(372)
D_{JK}	$\times 10^8$	3.26	3.23876(301)
H_J	$\times 10^{12}$	0.0	0.1846(343)
No. of data			729
σ (fit)	$\times 10^5$		5.1

^a Standard uncertainties (1σ) in parentheses refer to the least significant digits.

^b Ref. 24.

Overall, the metrological approach presented here considerably expands the portfolio of laser tools to probe molecules in the gas phase with high resolution and accuracy. This is proactive for the study of large molecules through the unique fingerprint features offered by the LW-MIR spectral region, for example molecules of the BTEX family (benzene, toluene, ethylbenzene, xylenes) that are of high interest for environmental modelling applications. In the case of benzene, which lacks a permanent dipole moment and, therefore, cannot be detected by pure rotational spectroscopy, infrared bands such as the ν_{11} explored here represent the clue to assess their presence in remote spatial environments, such as planetary atmospheres^{6,27,28}. On another front, by the use of a sub-hertz linewidth comb and phase-locking loops for pump and signal lasers, the spectrometer could reach the level of stability and spectral purity needed for tests of fundamental physics on molecular samples^{29,30} in a spectral region not explored so far.

METHODS

Difference Frequency Generation

The pump radiation in the DFG process is provided by a distributed-feedback quantum-cascade-laser (DFB-QCL, from Alpes Lasers) emitting around 5.69 μm . Tuning its temperature from 20 °C to -10 °C results in a continuous spectral coverage from 1757 to 1761 cm^{-1} (i.e., over 4 cm^{-1}) with an optical power from 17 to 48 mW, respectively. The signal radiation is generated by a CO₂ laser (L20GD, Access Laser) operating in a pulse modulation mode (“super-pulse mode”) with a duty cycle of 20 % and a modulation frequency of 10 kHz. Its emission can be finely piezo-tuned over few tens of MHz and coarsely tuned by grating over ~ 80 lines, from 9.23 to 10.71 μm , with peak powers varying from 3 to 27 W depending on the line. The wavelength combination of the two lasers provides access to 12.1-14.9 μm wavelength range for the idler beam, apart from 10 cm^{-1} gaps at 12.6 and 14.4 μm and a 30 cm^{-1} gap at 13.3 μm due to the discontinuities between adjacent CO₂ lasing bands. The OP-GaAs crystal used for DFG is 35 mm long and equipped with poling periods from 183 to 203 μm in a fan-out structure, which ensures phase matching over the entire tuning range available for the nonlinear laser source (upper bound at 14.8 μm rather than 14.9 μm). The width of the crystal is 20 mm, which is sufficiently large to discard any efficiency drop due to non-uniform phase-matching conditions across the beams. In fact, the poling period changes by less than 0.2 μm , which is far below the phase matching bandwidth of 0.7 μm , over beam diameters of 120 and 200 μm for pump and signal, respectively. As expressed in terms of wavenumbers, the phase-matching bandwidth amounts to 4.2 cm^{-1} , which implies that there is no need for adjustment of the crystal position while tuning the QCL. The optical power generated by the DFG considerably depends on both the QCL temperature and the CO₂ lasing line, but it remains above 20 μW (peak power) if the CO₂ emission is chosen among the 70 % most intense lines, as is the case for all measurements shown here.

Comb referencing and spectra acquisition

Pump (5.7 μm) and signal (9.2-10.7 μm) lasers are referenced to the same 1.9 μm Tm: fiber frequency comb through two independent sum frequency generation (SFG) processes, in AgGaSe₂ and ZnGaP (ZGP) crystals, respectively. The SFG processes shift the comb frequency $\nu_n = f_{\text{ceo}} + n f_{\text{rep}}$ (f_{ceo} and f_{rep} being carrier-envelope-offset and repetition frequency of the comb, respectively) by the pump (ν_p) and signal (ν_s) frequency, thereby generating two replicas of the comb around 1.4 ($\nu_{n,p} = \nu_n + \nu_p$) and 1.6 μm ($\nu_{n,s} = \nu_n + \nu_s$) for pump and signal, respectively. When these replicas are superimposed with a coherent continuum ($\nu_m = f_{\text{ceo}} + m f_{\text{rep}}$) generated from the original comb and extending down to 1.4 μm , a beat note is extracted, $f_{\text{beat}} = \pm(\nu_{n,p/s} - \nu_m) = \pm(\nu_{s,p} - (m-n)f_{\text{rep}})$ that straightforwardly links the unknown pump/signal frequency to an integer number of f_{rep} , independently of f_{ceo} . In our case, f_{rep} is stabilized to a GPS-tracked low-noise radiofrequency (RF) signal at 100 MHz while f_{ceo} is left free running. From an analysis of their beat notes

(Supplementary Fig. 1 and Fig. 2), QCL and CO₂ lasers give an almost Gaussian instrumental line shape of 1.6 and 1.5 MHz, respectively, which sums up to 2.2 MHz. The relatively large contribution from the CO₂ laser comes from the adopted super-pulse mode and it could be suppressed by choosing a cw regime followed by external modulation, at the expense of DFG power. The CO₂ beat note is locked to an RF local oscillator via piezo-feedback to the grating position, whereas the QCL beat note is measured in real time by a 100 MS/s DAQ, followed by an FFT conversion at every 1024 samples and a barycentre calculation via FPGA (PXIe-7961 FPGA board and NI-5781 add-on, National Instruments). The same board digitizes the lock-in output in order to synchronize horizontal and vertical axis of the measurement. An example of raw data for QCL beat note and gas absorption is given in the Supplementary Fig. 3. The acquisition of a 4 cm⁻¹ large spectrum is typically accomplished in 12 minutes (tuning rate of 180 MHz/s) to prevent laser frequency changes by more than its linewidth over the integration time of 10 ms set for the lock-in. A higher speed would be possible with a faster modulation of the CO₂ laser and a correspondingly smaller integration time. Spectra larger than 4 cm⁻¹, as those in Fig. 2 and Fig. 4, are acquired piecewise, due to the need of changing the CO₂ laser frequency and correspondingly adjusting the phase-matching condition for DFG. No tuning of the phase-matching angle is necessary for the SFG processes because of the much larger phase-matching bandwidths given by the shorter interaction length (6 mm) and the type of phase matching (birefringence rather than quasi-phase-matching). The assignment of the comb mode order for correct calibration of the frequency axis is straightforward because of the low uncertainty of both CO₂ and C₆H₆ lines, well below the comb mode spacing (100 MHz).

Uncertainty budget

A 240 MHz-large spectrum of a single isolated line, such as that shown in Fig. 2b, is typically acquired in 2 s and encompasses 170 spectral points: these are almost evenly distributed in 5 groups separated at every 50 MHz, namely at every half f_{rep} , by 18 MHz blanks due to the missed tracking of the comb-QCL beat note. In these conditions, the absorption noise ($3 \cdot 10^{-3}$ on the most intense lines) and the frequency uncertainty on each spectral point (29 kHz) account for about 70% and 5%, respectively, of the observed 140 kHz root mean square (rms) deviation of the line centre for repeated measurements. The missing contribution comes from distortions of the spectral baseline that change from measurement to measurement and are not properly fit out. The statistical uncertainty dominates over systematic effects, due to the negligible impact from the instability of the GPS-based frequency standard (2 kHz over 2 s), the inaccuracy of the absolute pressure gauge and of the pressure shift coefficients used to extrapolate line centres to zero pressure (< 2 kHz, thanks to the low pressure), and the inappropriateness of the Voigt fitting model (< 3 kHz, due to the negligible asymmetry of line profiles in the low-pressure collisional regime).

REFERENCES

1. Picqué, N. & Hänsch, T.W. Frequency comb spectroscopy. *Nat. Photonics* **146**, 146-157 (2019)
2. Baumann, E. & Fortier, T. 20 years of developments in optical frequency comb technology and applications. *Commun. Phys.* **2**, 153 (2019).
3. Jung, D., Bank, S., Lee, M. L. & Wasserman, D. Next-generation mid-infrared sources. *J. Opt.* **19**, 123001 (2017).
4. Argence, B et al. Quantum cascade laser frequency stabilisation at the sub-Hz level. *Nat. Photonics* **9**, 456-460 (2015).
5. Spaun, B. et al. Continuous probing of cold complex molecules with infrared frequency comb spectroscopy. *Nature* **533**, 517-520 (2016).
6. Vinatier, S. et al. Study of Titan's fall southern stratospheric polar cloud composition with Cassini/CIRS: detection of benzene ice. *Icarus* **310**, 89-104 (2018).

7. Consolino, L., Cappelli, F., de Cumis, M. S. & De Natale, P. QCL-based frequency metrology from the mid-infrared to the THz range: a review. *Nanophotonics* **8**, 181–204 (2019).
8. Tacke, M. Lead-salt lasers. *Philos. Trans. R. Soc. A: Mathematical, Physical and Engineering Sciences* **359**, 547–566 (2001).
9. Yao, Y., Hoffman, A.J. & Gmachl C.F. Mid-infrared quantum cascade lasers. *Nat. Photonics* **6**, 432–439 (2012).
10. Vitiello, M.S., Scalfari, G., Williams, B. & De Natale P. Quantum cascade lasers: 20 years of challenges. *Opt. Express* **23**, 5167–5182 (2015).
11. Nguyen Van, H. et al. Long Wavelength ($\lambda > 17 \mu\text{m}$) Distributed Feedback Quantum Cascade Lasers Operating in a Continuous Wave at Room Temperature. *Photonics*, **6**, 31 (2019).
12. Eckhoff, W., Putnam, R., Wang, S., Curl, R: & Tittel, F. A continuously tunable long-wavelength cw IR source for high resolution spectroscopy and trace gas detection. *Appl. Phys. B* **63**, 437–441 (1996).
13. Chen, W., Mouret, G. & Boucher, D. Difference-frequency laser spectroscopy detection of acetylene trace constituent. *Appl. Phys. B* **67**, 375–378 (1998).
14. Putnam, R.S., & Lancaster, D. G. Continuous-wave laser spectrometer automatically aligned and continuously tuned from 11.8 to 16.1 μm by use of diode-laser-pumped difference frequency generation in GaSe. *Appl. Opt.* **38**, 1513–1522 (1999).
15. Albert, S., Keppler, K., Quack, M. & Wokaun, A. Synchrotron-based highest resolution FTIR spectroscopy of chlorobenzene. *J. Mol. Spectrosc.* **315**, 92–101 (2015).
16. Cossel, K.C. et al. Gas-phase broadband spectroscopy using active sources: progress, status, and applications. *J. Opt. Soc. Am. B* **34**, 104–129 (2017).
17. Coddington, I., Newbury, N. & Swann, W. Dual-comb spectroscopy. *Optica* **3**, 414–426 (2016).
18. Kowligy, A.S. et al. Infrared electric field sampled frequency comb spectroscopy. *Sci. Adv.* **5** (2019).
19. Hansen, M.G., Magoulakis, E., Chen, Q.-F. & Schiller I.E.S. Quantum cascade laser-based mid-IR frequency metrology system with ultra-narrow linewidth and 1×10^{-13} -level frequency instability. *Opt. Lett.* **40**, 2289–2292 (2015).
20. Utsav K. C., Nasir, E.F. & Farooq, A. A mid-infrared absorption diagnostic for acetylene detection. *Appl. Phys. B* **120**, 223–232 (2015).
21. Schunemann, P.G. et al. Advances in nonlinear optical crystals for mid-infrared coherent sources. *J. Opt. Soc. Am. B* **33**, D36–D43 (2016).
22. Gambetta, A. et. al. Mid-infrared quantitative spectroscopy by comb-referencing of a quantum-cascade-laser: Application to the CO₂ spectrum at 4.3 μm . *Appl. Phys. Lett.* **99**, 251107 (2011).
23. Gordon, I. E. et al. The HITRAN2016 molecular spectroscopic database. *J. Quant. Spectrosc. Radiat. Transf.* **203**, 3–69 (2017).
24. Doi, A., Kasahara, S., Katô, & Baba M. Sub-Doppler rotationally resolved spectroscopy of lower vibronic bands of benzene with Zeeman effects. *J. Chem. Phys.* **120**, 6439–6448 (2004).
25. Wilson, E.B. The Statistical Weights of the Rotational Levels of Polyatomic Molecules, Including Methane, Ammonia, Benzene, Cyclopropane and Ethylene. *J. Chem. Phys.* **3**, 276–285 (1935).
26. Hollenstein, H., Piccirillo, S., Quack, M. & Snels, M. High-resolution infrared spectrum and analysis of the ν_{11} , A_{2u}(B₂) fundamental band of ¹²C₆H₆ and ¹³C¹²C₅H₆. *Mol. Phys.* **71**, 759–768 (1990).
27. Koskinen T.T., Moses J.I., West, R.A., Guerlet, S. & Jouchoux, A. The detection of benzene in Saturn’s upper atmosphere. *Geophysical Res. Lett.* **43**, 7895–7901 (2016).
28. Bézard, B., Drossart, P. & Encrenaz T. Benzene on the giant planets. *Icarus* **154**, 492–500 (2001).
29. Baron, J. et al. Order of magnitude smaller limit on the electric dipole moment of the electron. *Science* **343**, 269–272 (2014).
30. Daussy, C. et al. First direct determination of the Boltzmann constant by an optical method: towards a new definition of the kelvin. *Phys. Rev. Lett.* **98**, 250801 (2007).

ACKNOWLEDGMENTS

The authors acknowledge a financial contribution from the cooperative project OSR-2019-CCF-1975.34 between Politecnico di Milano and King Abdullah University of Science and Technology (KAUST) and by the project EMPATIA@Lecco ID: 2016-1428. F.T. and E.C. thank the University of Bologna for RFO funds.

AUTHOR CONTRIBUTIONS

M.M. and A.F. conceived the experiments. M.L and K.S. realized and characterized the DFG source while D.G. and R.G. took care of its referencing to the comb. M.L., R.G., D.G. performed the spectroscopic measurements. F.T., E.C. took care of fitting and interpreting benzene spectra, while M.L., R.G. and P.L. of carbon dioxide spectra. P.S. designed and provided the nonlinear crystal. M.M. took care of the first draft writing. All authors contributed to and edited the manuscript.

SUPPLEMENTARY MATERIAL

Bending modes metrology in the 14-15 μm region

M. Lamperti¹, R. Gotti¹, D. Gatti¹, M. K. Shakfa², E. Cané³, F. Tamassia³, P. Schunemann⁴, P. Laporta¹, A. Farooq^{2*}, and M. Marangoni^{1*}

¹Dipartimento di Fisica - Politecnico di Milano and IFN-CNR, Via Gaetano Prevati 1/C, 23900 Lecco, Italy

²King Abdullah University for Science and Technology, Clean Combustion Research Center, Thuwal 23955, Saudi Arabia

³Università di Bologna, Dipartimento di Chimica Industriale, Viale del Risorgimento 4, 40136 Bologna, Italy

⁴BAE Systems, Inc., MER15-1813, P.O. Box 868, Nashua, New Hampshire 03061-0868, USA

*amir.farooq@kaust.edu.sa, marco.marangoni@polimi.it

Characterization of the DFG linewidth

The linewidth of the DFG laser source has been estimated from the measured beat notes of the QCL and CO₂ lasers with the comb (Fig. S1 and S2, respectively). These have an almost Gaussian-like shape, which justifies a quadrature sum of their widths to have a reliable estimation of the DFG linewidth. The contribution of the comb linewidth (~ 100 kHz) to the beat notes (~ 1.5 MHz) is negligible.

For both lasers, the beat note has been measured via the 100 MS/s DAQ used for spectral acquisition, by subdividing the beat note samples in groups of 1024 points and calculating the barycentre of the electrical spectrum for each group via FFT. This is equivalent to repeatedly measuring the average laser frequency (as compared to the comb) over a measurement time of 10 μs (1024 points at every 10 ns). The distribution of these frequencies is shown in Fig. S1 and S2 for QCL and CO₂ laser, respectively. The width of these distributions provides an estimation value for their linewidths.

Supplementary Fig. 1

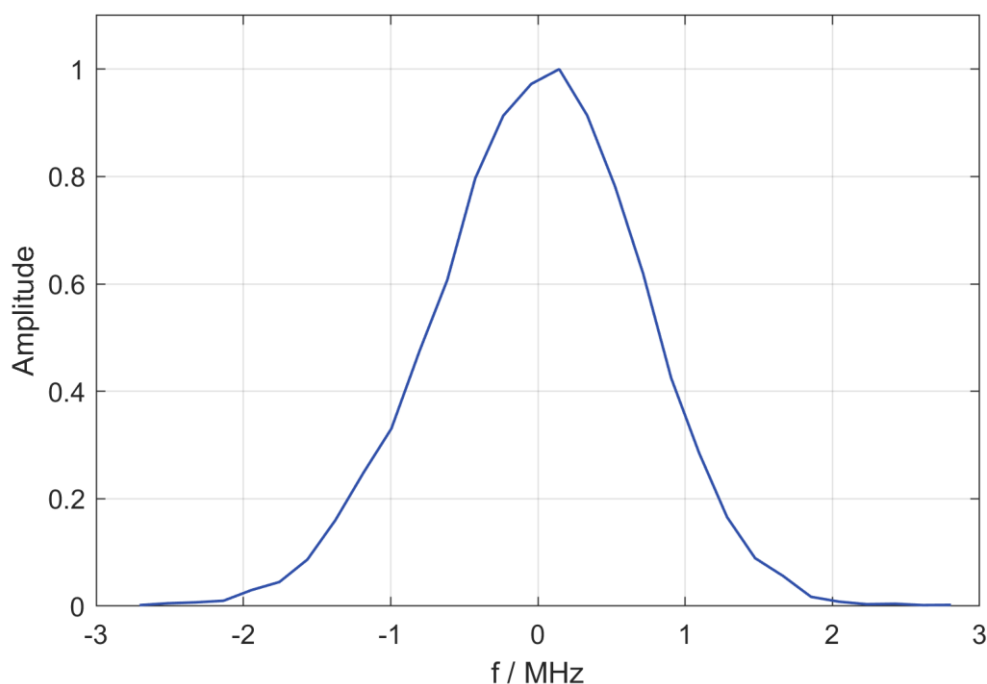


Fig. S1 | Beat note spectrum between QCL and comb. It is obtained as a histogram of the barycentre frequencies measured by the acquisition board over subsequent time windows of 10 μs .

Supplementary Fig. 2

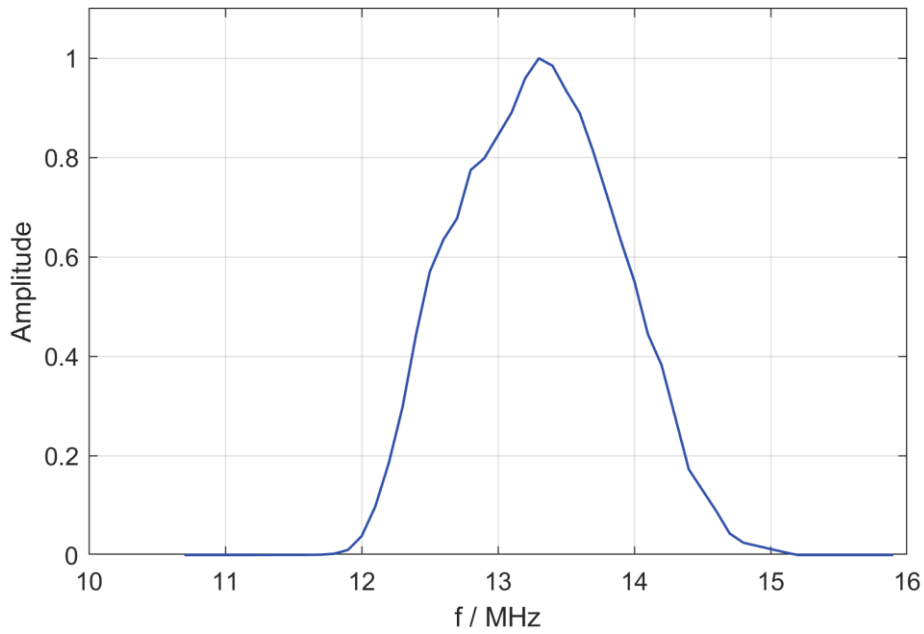


Fig. S2 | Beat note spectrum between CO₂ laser and comb. It is obtained as a histogram of the barycentre frequencies measured by the acquisition board over subsequent time windows of 10 μ s. The slight asymmetry derives from the pulsed operation of the laser.

Supplementary Fig. 3

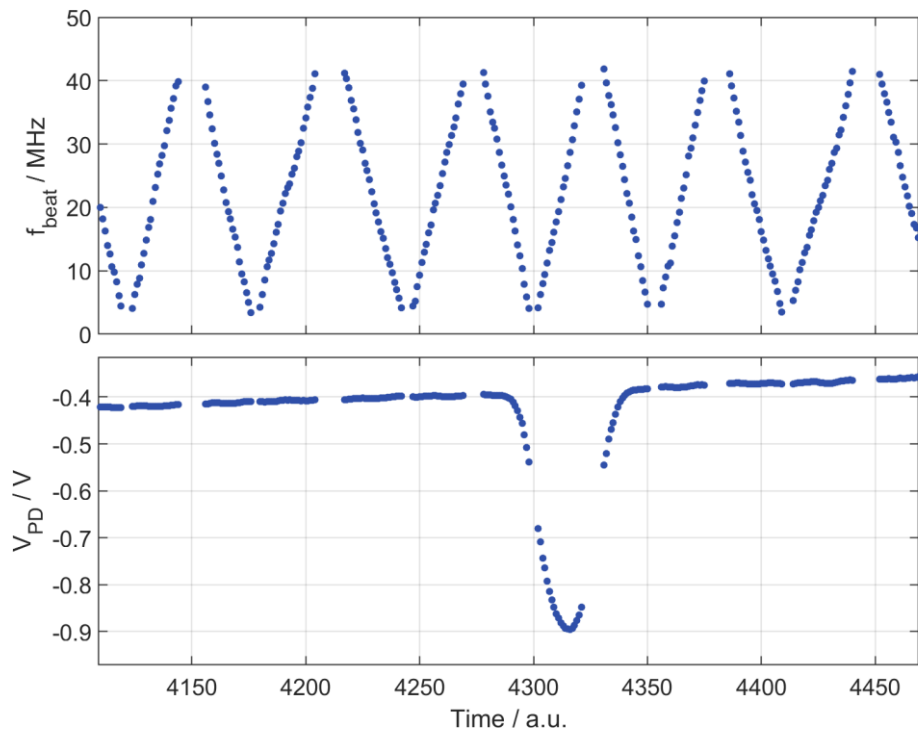


Fig. S3 | Raw measurement data for a CO₂ line. Time behaviour of the beat note frequency between QCL and comb (upper panel) and of the photodetector signal (lower panel).

Supplementary Table 1**Table S1 | Carbon dioxide line list.**

Observed and calculated wavenumbers of the transitions assigned to ν_2 and $\nu_1 \leftarrow \nu_2$ for CO_2 , ordered for increasing energy.

Band	Line	Obs / cm^{-1}	HITRAN/ cm^{-1}	Obs-HITRAN / MHz
ν_2	R(10)	676.01943603(66)	676.0194310	0.151
ν_2	R(12)	677.6008247(13)	677.6008210	0.113
ν_2	R(14)	679.18538721(70)	679.1853870	0.006
ν_2	R(16)	680.7731066(15)	680.7730990	0.228
ν_2	R(18)	682.3639287(18)	682.3639270	0.052
ν_2	R(20)	683.9578497(10)	683.9578360	0.413
ν_2	R(22)	685.5548060(34)	685.5547940	0.360
ν_2	R(24)	687.1547798(30)	687.1547670	0.386
ν_2	R(26)	688.7577330(62)	688.7577190	0.422
$\nu_1 \leftarrow \nu_2$	Q(36)	719.4226506(45)	719.4226250	0.769
$\nu_1 \leftarrow \nu_2$	Q(34)	719.5662865(32)	719.5662710	0.466
$\nu_1 \leftarrow \nu_2$	Q(32)	719.7026444(29)	719.7026280	0.494
$\nu_1 \leftarrow \nu_2$	Q(30)	719.8315747(41)	719.8315570	0.532
$\nu_1 \leftarrow \nu_2$	Q(28)	719.9529553(39)	719.9529220	1.000
$\nu_1 \leftarrow \nu_2$	Q(26)	720.0666311(55)	720.0666010	0.905
$\nu_1 \leftarrow \nu_2$	Q(24)	720.1724819(58)	720.1724780	0.118
$\nu_1 \leftarrow \nu_2$	Q(22)	720.2704506(72)	720.2704460	0.139
$\nu_1 \leftarrow \nu_2$	Q(18)	720.4422813(75)	720.4422760	0.160
$\nu_1 \leftarrow \nu_2$	Q(14)	720.5814236(50)	720.5814110	0.378
$\nu_1 \leftarrow \nu_2$	Q(12)	720.6385565(71)	720.6385440	0.377
$\nu_1 \leftarrow \nu_2$	Q(8)	720.7276548(51)	720.7276640	-0.276
$\nu_1 \leftarrow \nu_2$	Q(6)	720.7595686(35)	720.7595650	0.110
$\nu_1 \leftarrow \nu_2$	Q(4)	720.7829825(26)	720.7829840	-0.045
$\nu_1 \leftarrow \nu_2$	Q(2)	720.7979103(34)	720.7978970	0.401

Supplementary Table 2

Table S2 | Benzene line list.

Observed and calculated wavenumbers of the transitions assigned to ν_{11} for benzene, ordered for increasing energy.

The headings J', K' and J'', K'' indicate the values of the rotational quantum numbers for the upper and lower states of transitions, respectively. The weights used in the fitting of band constants are also reported (see the main text for description).

J'	K'	J''	K''	OBS / cm ⁻¹	CALC / cm ⁻¹	OBS-CALC / cm ⁻¹	WEIGHT
3	0	2	0	675.11231800	675.11215451	0.00016349	0.0033
3	1	2	1	675.11231800	675.11232569	-0.00000769	0.0033
3	2	2	2	675.11231800	675.11283922	-0.00052122	0.0033
4	0	3	0	675.49087000	675.49061191	0.00025809	0.0033
4	2	3	2	675.49087000	675.49129713	-0.00042713	0.0033
4	1	3	1	675.49087000	675.49078321	0.00008679	0.0033
4	3	3	3	675.49209300	675.49215368	-0.00006068	0.0100
5	0	4	0	675.86904600	675.86879520	0.00025080	0.0033
5	2	4	2	675.86904600	675.86948094	-0.00043494	0.0033
5	1	4	1	675.86904600	675.86896663	0.00007937	0.0033
5	3	4	3	675.87012900	675.87033813	-0.00020913	0.0100
5	4	4	4	675.87141300	675.87153821	-0.00012521	0.0100
6	0	5	0	676.24695800	676.24670350	0.00025450	0.0033
6	2	5	2	676.24695800	676.24738976	-0.00043176	0.0033
6	1	5	1	676.24695800	676.24687506	0.00008294	0.0033
6	3	5	3	676.24807800	676.24824759	-0.00016959	0.0100
6	4	5	4	676.24943000	676.24944857	-0.00001857	0.0100
6	5	5	5	676.25094400	676.25099272	-0.00004872	0.0100
7	0	6	0	676.62457400	676.62433591	0.00023809	0.0033
7	1	6	1	676.62457400	676.62450760	0.00006640	0.0033
7	2	6	2	676.62457400	676.62502268	-0.00044868	0.0033
7	3	6	3	676.62567700	676.62588115	-0.00020415	0.0100
7	4	6	4	676.62702000	676.62708303	-0.00006303	0.0100
7	5	6	5	676.62870788	676.62862833	0.00007955	1.0000

7	6	6	6	676.63046360	676.63051708	-0.00005348	1.0000
8	0	7	0	677.00197100	677.00169150	0.00027950	0.0033
8	1	7	1	677.00197100	677.00186332	0.00010768	0.0033
8	2	7	2	677.00197100	677.00237878	-0.00040778	0.0033
8	3	7	3	677.00310700	677.00323789	-0.00013089	0.0100
8	4	7	4	677.00426100	677.00444066	-0.00017966	0.0100
8	5	7	5	677.00599574	677.00598711	0.00000863	1.0000
8	6	7	6	677.00791100	677.00787727	0.00003373	0.0100
8	7	7	7	677.01003686	677.01011116	-0.00007430	1.0000
9	0	8	0	677.37899200	677.37876935	0.00022265	0.0033
9	1	8	1	677.37899200	677.37894129	0.00005071	0.0033
9	2	8	2	677.37899200	677.37945714	-0.00046514	0.0033
9	3	8	3	677.38015700	677.38031688	-0.00015988	0.0100
9	4	8	4	677.38147600	677.38152054	-0.00004454	0.0100
9	5	8	5	677.38308347	677.38306814	0.00001533	1.0000
9	6	8	6	677.38495406	677.38495970	-0.00000564	1.0000
9	7	8	7	677.38717956	677.38719523	-0.00001567	1.0000
9	8	8	8	677.38975483	677.38977479	-0.00001996	1.0000
10	0	9	0	677.75584400	677.75556850	0.00027550	0.0033
10	2	9	2	677.75584400	677.75625680	-0.00041280	0.0033
10	1	9	1	677.75584400	677.75574058	0.00010342	0.0033
10	3	9	3	677.75703300	677.75711718	-0.00008418	0.0100
10	4	9	4	677.75812100	677.75832173	-0.00020073	0.0100
10	5	9	5	677.75987231	677.75987047	0.00000184	1.0000
10	6	9	6	677.76177450	677.76176341	0.00001109	1.0000
10	7	9	7	677.76396400	677.76400060	-0.00003660	0.0100
10	8	9	8	677.76659887	677.76658206	0.00001681	1.0000
10	9	9	9	677.76950159	677.76950783	-0.00000624	1.0000
11	0	10	0	678.13232400	678.13208801	0.00023599	0.0033
11	2	10	2	678.13232400	678.13277681	-0.00045281	0.0033
11	1	10	1	678.13232400	678.13226021	0.00006379	0.0033

11	3	10	3	678.13340500	678.13363782	-0.00023282	0.0100
11	4	10	4	678.13481400	678.13484326	-0.00002926	0.0100
11	5	10	5	678.13640025	678.13639313	0.00000712	1.0000
11	6	10	6	678.13826997	678.13828747	-0.00001750	1.0000
11	7	10	7	678.14051568	678.14052629	-0.00001061	1.0000
11	8	10	8	678.14310128	678.14310965	-0.00000837	1.0000
11	9	10	9	678.14602413	678.14603756	-0.00001343	1.0000
11	10	10	10	678.14934367	678.14931008	0.00003359	1.0000
12	0	11	0	678.50862300	678.50832691	0.00029609	0.0033
12	2	11	2	678.50862300	678.50901621	-0.00039321	0.0033
12	1	11	1	678.50862300	678.50849923	0.00012377	0.0033
12	3	11	3	678.50975100	678.50987785	-0.00012685	0.0100
12	4	11	4	678.51105500	678.51108416	-0.00002916	0.0100
12	5	11	5	678.51265241	678.51263517	0.00001724	1.0000
12	6	11	6	678.51452983	678.51453089	-0.00000106	1.0000
12	7	11	7	678.51674900	678.51677135	-0.00002235	0.0100
12	8	11	8	678.51932347	678.51935659	-0.00003312	1.0000
12	9	11	9	678.52226729	678.52228664	-0.00001935	1.0000
12	10	11	10	678.52551193	678.52556155	-0.00004962	1.0000
12	11	11	11	678.52927988	678.52918136	0.00009852	1.0000
13	0	12	0	678.88451800	678.88428421	0.00023379	0.0033
13	2	12	2	678.88451800	678.88497401	-0.00045601	0.0033
13	1	12	1	678.88451800	678.88445666	0.00006134	0.0033
13	3	12	3	678.88558400	678.88583628	-0.00025228	0.0100
13	4	12	4	678.88699400	678.88704347	-0.00004947	0.0100
13	5	12	5	678.88860537	678.88859560	0.00000977	1.0000
13	6	12	6	678.89049281	678.89049270	0.00000011	1.0000
13	7	12	7	678.89274651	678.89273479	0.00001172	1.0000
13	8	12	8	678.89529787	678.89532191	-0.00002404	1.0000
13	9	12	9	678.89824600	678.89825409	-0.00000809	0.0100
13	10	12	10	678.90151000	678.90153138	-0.00002138	0.0100

13	11	12	11	678.90513200	678.90515382	-0.00002182	0.0100
13	12	12	12	678.91000700	678.90912148	0.00088552	0.0050
14	0	13	0	679.26030500	679.25995893	0.00034607	0.0033
14	2	13	2	679.26030500	679.26064923	-0.00034423	0.0033
14	1	13	1	679.26030500	679.26013151	0.00017349	0.0033
14	3	13	3	679.26140700	679.26151212	-0.00010512	0.0100
14	4	13	4	679.26257100	679.26272019	-0.00014919	0.0100
14	5	13	5	679.26428299	679.26427344	0.00000955	1.0000
14	6	13	6	679.26616363	679.26617192	-0.00000829	1.0000
14	7	13	7	679.26840700	679.26841563	-0.00000863	1.0000
14	8	13	8	679.27100571	679.27100462	0.00000109	1.0000
14	9	13	9	679.27393347	679.27393892	-0.00000545	1.0000
14	10	13	10	679.27720541	679.27721858	-0.00001317	1.0000
14	11	13	11	679.28083430	679.28084365	-0.00000935	1.0000
14	12	13	12	679.28489600	679.28481418	0.00008182	0.0100
14	13	13	13	679.28925900	679.28913022	0.00012878	0.0100
15	0	14	0	679.63559600	679.63535009	0.00024591	0.0033
15	2	14	2	679.63559600	679.63604089	-0.00044489	0.0033
15	1	14	1	679.63559600	679.63552279	0.00007321	0.0033
15	3	14	3	679.63680700	679.63690440	-0.00009740	0.0100
15	4	14	4	679.63797000	679.63811334	-0.00014334	0.0100
15	5	14	5	679.63967371	679.63966771	0.00000600	1.0000
15	6	14	6	679.64157289	679.64156755	0.00000534	1.0000
15	7	14	7	679.64380547	679.64381288	-0.00000741	1.0000
15	8	14	8	679.64639884	679.64640373	-0.00000489	1.0000
15	9	14	9	679.64932588	679.64934015	-0.00001427	1.0000
15	10	14	10	679.65260752	679.65262218	-0.00001466	1.0000
15	11	14	11	679.65623920	679.65624985	-0.00001065	1.0000
15	12	14	12	679.65976600	679.66022324	-0.00045724	0.0100
15	13	14	13	679.66427800	679.66454239	-0.00026439	0.0100
15	14	14	14	679.66901700	679.66920738	-0.00019038	0.0100

16	0	15	0	680.01078500	680.01045669	0.00032831	0.0033
16	1	15	1	680.01078500	680.01062951	0.00015549	0.0033
16	2	15	2	680.01078500	680.01114798	-0.00036298	0.0033
16	3	15	3	680.01190000	680.01201211	-0.00011211	0.0100
16	4	15	4	680.01307900	680.01322192	-0.00014292	0.0100
16	5	15	5	680.01479009	680.01477741	0.00001268	1.0000
16	6	15	6	680.01669515	680.01667861	0.00001654	1.0000
16	7	15	7	680.01892957	680.01892555	0.00000402	1.0000
16	8	15	8	680.02149780	680.02151826	-0.00002046	1.0000
16	9	15	9	680.02445884	680.02445678	0.00000206	1.0000
16	10	15	10	680.02774118	680.02774116	0.00000002	1.0000
16	11	15	11	680.03133613	680.03137144	-0.00003531	1.0000
16	12	15	12	680.03536700	680.03534768	0.00001932	0.0100
16	13	15	13	680.03972848	680.03966993	0.00005855	1.0000
16	14	15	14	680.04446039	680.04433826	0.00012213	1.0000
16	15	15	15	680.04920015	680.04935273	-0.00015258	1.0000
17	0	16	0	680.38557700	680.38527773	0.00029927	0.0033
17	1	16	1	680.38557700	680.38545068	0.00012632	0.0033
17	2	16	2	680.38557700	680.38596952	-0.00039252	0.0033
17	3	16	3	680.38661700	680.38683427	-0.00021727	0.0100
17	4	16	4	680.38797700	680.38804493	-0.00006793	0.0100
17	5	16	5	680.38960629	680.38960153	0.00000476	1.0000
17	6	16	6	680.39149895	680.39150409	-0.00000514	1.0000
17	7	16	7	680.39374752	680.39375264	-0.00000512	1.0000
17	8	16	8	680.39635941	680.39634720	0.00001221	1.0000
17	9	16	9	680.39928742	680.39928782	-0.00000040	1.0000
17	10	16	10	680.40257273	680.40257455	-0.00000182	1.0000
17	11	16	11	680.40617662	680.40620742	-0.00003080	1.0000
17	12	16	12	680.41019317	680.41018649	0.00000668	1.0000
17	13	16	13	680.41452322	680.41451183	0.00001139	1.0000
17	14	16	14	680.41922800	680.41918349	0.00004451	0.0100

17	15	16	15	680.42423515	680.42420154	0.00003361	1.0000
17	16	16	16	680.42981802	680.42956606	0.00025196	1.0000
18	0	17	0	680.76015100	680.75981222	0.00033878	0.0033
18	1	17	1	680.76015100	680.75998529	0.00016571	0.0033
18	2	17	2	680.76015100	680.76050450	-0.00035350	0.0033
18	3	17	3	680.76127800	680.76136987	-0.00009187	0.0100
18	4	17	4	680.76257100	680.76258139	-0.00001039	0.0100
18	5	17	5	680.76416288	680.76413910	0.00002378	1.0000
18	6	17	6	680.76604703	680.76604301	0.00000402	1.0000
18	7	17	7	680.76828530	680.76829315	-0.00000785	1.0000
18	8	17	8	680.77089814	680.77088956	0.00000858	1.0000
18	9	17	9	680.77383052	680.77383227	-0.00000175	1.0000
18	10	17	10	680.77712629	680.77712133	0.00000496	1.0000
18	11	17	11	680.78098300	680.78075679	0.00022621	0.0100
18	12	17	12	680.78474486	680.78473869	0.00000617	1.0000
18	13	17	13	680.78900540	680.78906710	-0.00006170	1.0000
18	14	17	14	680.79370661	680.79374208	-0.00003547	1.0000
18	15	17	15	680.79874459	680.79876370	-0.00001911	1.0000
18	16	17	16	680.80396663	680.80413203	-0.00016540	1.0000
18	17	17	17	680.81001348	680.80984715	0.00016633	1.0000
19	0	18	0	681.13434200	681.13405917	0.00028283	0.0033
19	2	18	2	681.13434200	681.13475194	-0.00040994	0.0033
19	1	18	1	681.13434200	681.13423236	0.00010964	0.0033
19	3	18	3	681.13543700	681.13561792	-0.00018092	0.0100
19	4	18	4	681.13675900	681.13683030	-0.00007130	0.0100
19	5	18	5	681.13839905	681.13838911	0.00000994	1.0000
19	6	18	6	681.14028400	681.14029437	-0.00001037	0.0100
19	7	18	7	681.14253884	681.14254611	-0.00000727	1.0000
19	8	18	8	681.14512170	681.14514435	-0.00002265	1.0000
19	9	18	9	681.14810284	681.14808914	0.00001370	1.0000
19	10	18	10	681.15137798	681.15138053	-0.00000255	1.0000

19	11	18	11	681.15494733	681.15501855	-0.00007122	1.0000
19	12	18	12	681.15907625	681.15900327	0.00007298	1.0000
19	13	18	13	681.16336997	681.16333474	0.00003523	1.0000
19	14	18	14	681.16800786	681.16801303	-0.00000517	1.0000
19	15	18	15	681.17310329	681.17303820	0.00006509	1.0000
19	16	18	16	681.17839418	681.17841033	-0.00001615	1.0000
19	17	18	17	681.18403614	681.18412950	-0.00009336	1.0000
19	18	18	18	681.19023828	681.19019579	0.00004249	1.0000
20	0	19	0	681.50836800	681.50801759	0.00035041	0.0033
20	2	19	2	681.50836800	681.50871085	-0.00034285	0.0033
20	1	19	1	681.50836800	681.50819091	0.00017709	0.0033
20	3	19	3	681.50954200	681.50957744	-0.00003544	0.0100
20	4	19	4	681.51070200	681.51079067	-0.00008867	0.0100
20	5	19	5	681.51239692	681.51235058	0.00004634	1.0000
20	6	19	6	681.51427623	681.51425718	0.00001905	1.0000
20	7	19	7	681.51652533	681.51651050	0.00001483	1.0000
20	8	19	8	681.51909555	681.51911058	-0.00001503	1.0000
20	9	19	9	681.52204404	681.52205745	-0.00001341	1.0000
20	10	19	10	681.52534503	681.52535115	-0.00000612	1.0000
20	11	19	11	681.52893300	681.52899174	-0.00005874	0.0050
20	12	19	12	681.53295413	681.53297926	-0.00002513	1.0000
20	13	19	13	681.53730080	681.53731378	-0.00001298	1.0000
20	14	19	14	681.54200654	681.54199536	0.00001118	1.0000
20	15	19	15	681.54697987	681.54702407	-0.00004420	1.0000
20	16	19	16	681.55249291	681.55239998	0.00009293	1.0000
20	17	19	17	681.55823512	681.55812318	0.00011194	1.0000
20	18	19	18	681.56375842	681.56419374	-0.00043532	1.0000
20	19	19	19	681.57073200	681.57061175	0.00012025	0.0100
21	0	20	0	681.88188600	681.88168650	0.00019950	0.0033
21	2	20	2	681.88188600	681.88238025	-0.00049425	0.0033
21	1	20	1	681.88188600	681.88185994	0.00002606	0.0033

21	3	20	3	681.88310900	681.88324744	-0.00013844	0.0100
21	4	20	4	681.88437700	681.88446153	-0.00008453	0.0100
21	5	20	5	681.88607339	681.88602253	0.00005086	1.0000
21	6	20	6	681.88795351	681.88793047	0.00002304	1.0000
21	7	20	7	681.89018781	681.89018537	0.00000244	1.0000
21	8	20	8	681.89279003	681.89278727	0.00000276	1.0000
21	9	20	9	681.89574156	681.89573620	0.00000536	1.0000
21	10	20	10	681.89903791	681.89903221	0.00000570	1.0000
21	11	20	11	681.90253000	681.90267535	-0.00014535	0.0100
21	12	20	12	681.90669123	681.90666567	0.00002556	1.0000
21	13	20	13	681.91098391	681.91100323	-0.00001932	1.0000
21	14	20	14	681.91579495	681.91568809	0.00010686	1.0000
21	15	20	15	681.92067553	681.92072033	-0.00004480	1.0000
21	16	20	16	681.92601896	681.92610001	-0.00008105	1.0000
21	17	20	17	681.93179117	681.93182721	-0.00003604	1.0000
21	18	20	18	681.93791084	681.93790202	0.00000882	1.0000
21	19	20	19	681.94438026	681.94432453	0.00005573	1.0000
21	20	20	20	681.95107400	681.95109484	-0.00002084	0.0050
22	0	21	0	682.25544200	682.25506493	0.00037707	0.0033
22	1	21	1	682.25544200	682.25523849	0.00020351	0.0033
22	2	21	2	682.25544200	682.25575916	-0.00031716	0.0033
22	3	21	3	682.25647600	682.25662696	-0.00015096	0.0050
22	4	21	4	682.25784100	682.25784190	-0.00000090	0.0100
22	5	21	5	682.25943765	682.25940399	0.00003366	1.0000
22	6	21	6	682.26133265	682.26131326	0.00001939	1.0000
22	7	21	7	682.26359938	682.26356973	0.00002965	1.0000
22	8	21	8	682.26616458	682.26617345	-0.00000887	1.0000
22	9	21	9	682.26912962	682.26912444	0.00000518	1.0000
22	10	21	10	682.27248200	682.27242275	0.00005925	0.0100
22	11	21	11	682.27613390	682.27606843	0.00006547	1.0000
22	12	21	12	682.28021600	682.28006153	0.00015447	0.0100

22	13	21	13	682.28443414	682.28440212	0.00003202	1.0000
22	14	21	14	682.28905663	682.28909025	-0.00003362	1.0000
22	15	21	15	682.29416584	682.29412600	0.00003984	1.0000
22	16	21	16	682.29957031	682.29950943	0.00006088	1.0000
22	17	21	17	682.30528747	682.30524063	0.00004684	1.0000
22	18	21	18	682.31140752	682.31131968	0.00008784	1.0000
22	19	21	19	682.31762569	682.31774667	-0.00012098	1.0000
22	20	21	20	682.32468219	682.32452169	0.00016050	1.0000
22	21	21	21	682.33139871	682.33164485	-0.00024614	1.0000
23	0	22	0	682.62844200	682.62815193	0.00029007	0.0033
23	2	22	2	682.62844200	682.62884664	-0.00040464	0.0033
23	1	22	1	682.62844200	682.62832560	0.00011640	0.0033
23	3	22	3	682.62957100	682.62971504	-0.00014404	0.0100
23	4	22	4	682.63092600	682.63093082	-0.00000482	0.0100
23	5	22	5	682.63249998	682.63249400	0.00000598	1.0000
23	6	22	6	682.63440495	682.63440459	0.00000036	1.0000
23	7	22	7	682.63665849	682.63666263	-0.00000414	1.0000
23	8	22	8	682.63926055	682.63926816	-0.00000761	1.0000
23	10	22	10	682.64561771	682.64552180	0.00009591	1.0000
23	11	22	11	682.64919153	682.64917002	0.00002151	1.0000
23	12	22	12	682.65307320	682.65316589	-0.00009269	1.0000
23	13	22	13	682.65752704	682.65750949	0.00001755	1.0000
23	14	22	14	682.66222210	682.66220088	0.00002122	1.0000
23	15	22	15	682.66721275	682.66724012	-0.00002737	1.0000
23	16	22	16	682.67253200	682.67262729	-0.00009529	0.0100
23	17	22	17	682.67837215	682.67836247	0.00000968	1.0000
23	18	22	18	682.68446923	682.68444574	0.00002349	1.0000
23	19	22	19	682.69104700	682.69087719	0.00016981	0.0100
23	20	22	20	682.69760017	682.69765692	-0.00005675	1.0000
23	21	22	21	682.70484037	682.70478503	0.00005534	1.0000
23	22	22	22	682.71235300	682.71226161	0.00009139	0.0100

24	0	23	0	683.00135200	683.00094654	0.00040546	0.0033
24	1	23	1	683.00135200	683.00112034	0.00023166	0.0033
24	2	23	2	683.00135200	683.00164174	-0.00028974	0.0033
24	3	23	3	683.00247500	683.00251074	-0.00003574	0.0100
24	4	23	4	683.00371000	683.00372736	-0.00001736	0.0100
24	5	23	5	683.00531188	683.00529162	0.00002026	1.0000
24	6	23	6	683.00719380	683.00720353	-0.00000973	1.0000
24	7	23	7	683.00947277	683.00946314	0.00000963	1.0000
24	8	23	8	683.01206706	683.01207046	-0.00000340	1.0000
24	9	23	9	683.01500308	683.01502555	-0.00002247	1.0000
24	10	23	10	683.01818942	683.01832843	-0.00013901	1.0000
24	11	23	11	683.02205123	683.02197917	0.00007206	1.0000
24	12	23	12	683.02606217	683.02597781	0.00008436	1.0000
24	13	23	13	683.03030803	683.03032441	-0.00001638	1.0000
24	14	23	14	683.03503423	683.03501904	0.00001519	1.0000
24	15	23	15	683.04009338	683.04006177	0.00003161	1.0000
24	16	23	16	683.04544511	683.04545266	-0.00000755	1.0000
24	17	23	17	683.05121836	683.05119180	0.00002656	1.0000
24	18	23	18	683.05731608	683.05727928	0.00003680	1.0000
24	19	23	19	683.06374111	683.06371517	0.00002594	1.0000
24	20	23	20	683.07060100	683.07049959	0.00010141	0.0100
24	21	23	21	683.07765430	683.07763262	0.00002168	1.0000
24	22	23	22	683.08523321	683.08511436	0.00011885	1.0000
24	23	23	23	683.09298327	683.09294494	0.00003833	1.0000
25	0	24	0	683.37364600	683.37344786	0.00019814	0.0033
25	2	24	2	683.37364600	683.37414353	-0.00049753	0.0033
25	1	24	1	683.37364600	683.37362178	0.00002422	0.0033
25	3	24	3	683.37496400	683.37501313	-0.00004913	0.0100
25	4	24	4	683.37616300	683.37623059	-0.00006759	0.0100
25	5	24	5	683.37780764	683.37779592	0.00001172	1.0000
25	6	24	6	683.37970959	683.37970916	0.00000043	1.0000

25	7	24	7	683.38196299	683.38197032	-0.00000733	1.0000
25	8	24	8	683.38457566	683.38457944	-0.00000378	1.0000
25	9	24	9	683.38752977	683.38753655	-0.00000678	1.0000
25	10	24	10	683.39090800	683.39084171	0.00006629	0.0100
25	11	24	11	683.39449765	683.39449496	0.00000269	1.0000
25	12	24	12	683.39849753	683.39849635	0.00000118	1.0000
25	13	24	13	683.40285846	683.40284595	0.00001251	1.0000
25	14	24	14	683.40752791	683.40754381	-0.00001590	1.0000
25	15	24	15	683.41257082	683.41259000	-0.00001918	1.0000
25	16	24	16	683.41798711	683.41798461	0.00000250	1.0000
25	17	24	17	683.42369478	683.42372770	-0.00003292	1.0000
25	18	24	18	683.42988073	683.42981936	0.00006137	1.0000
25	19	24	19	683.43622350	683.43625969	-0.00003619	1.0000
25	20	24	20	683.44304319	683.44304876	-0.00000557	1.0000
25	21	24	21	683.45016853	683.45018670	-0.00001817	1.0000
25	22	24	22	683.45775731	683.45767359	0.00008372	1.0000
25	23	24	23	683.46547324	683.46550955	-0.00003631	1.0000
25	24	24	24	683.47380800	683.47369470	0.00011330	0.0100
26	0	25	0	683.74609400	683.74565498	0.00043902	0.0033
26	2	25	2	683.74609400	683.74635113	-0.00025713	0.0033
26	1	25	1	683.74609400	683.74582901	0.00026499	0.0033
26	3	25	3	683.74705600	683.74722132	-0.00016532	0.0100
26	4	25	4	683.74841300	683.74843961	-0.00002661	0.0100
26	5	25	5	683.75001819	683.75000602	0.00001217	1.0000
26	6	25	6	683.75191750	683.75192056	-0.00000306	1.0000
26	7	25	7	683.75416577	683.75418327	-0.00001750	1.0000
26	8	25	8	683.75679197	683.75679418	-0.00000221	1.0000
26	9	25	9	683.75973952	683.75975332	-0.00001380	1.0000
26	10	25	10	683.76304951	683.76306075	-0.00001124	1.0000
26	11	25	11	683.76668485	683.76671650	-0.00003165	1.0000
26	12	25	12	683.77070222	683.77072063	-0.00001841	1.0000

26	13	25	13	683.77510961	683.77507321	0.00003640	1.0000
26	14	25	14	683.77977672	683.77977429	0.00000243	1.0000
26	15	25	15	683.78483947	683.78482394	0.00001553	1.0000
26	16	25	16	683.79023298	683.79022224	0.00001074	1.0000
26	17	25	17	683.79606170	683.79596926	0.00009244	1.0000
26	18	25	18	683.80202731	683.80206509	-0.00003778	1.0000
26	19	25	19	683.80851600	683.80850983	0.00000617	1.0000
26	20	25	20	683.81537300	683.81530355	0.00006945	0.0100
26	21	25	21	683.82249035	683.82244637	0.00004398	1.0000
26	22	25	22	683.82988927	683.82993839	-0.00004912	1.0000
26	23	25	23	683.83787300	683.83777972	0.00009328	0.0100
26	24	25	24	683.84612298	683.84597046	0.00015252	1.0000
26	25	25	25	683.85425121	683.85451075	-0.00025954	1.0000
27	0	26	0	684.11774700	684.11756702	0.00017998	0.1667
27	2	26	2	684.11774700	684.11826364	-0.00051664	0.0017
27	1	26	1	684.11774700	684.11774118	0.00000582	0.0017
27	3	26	3	684.11905400	684.11913443	-0.00008043	0.0100
27	4	26	4	684.12022800	684.12035356	-0.00012556	0.0100
27	5	26	5	684.12191637	684.12192103	-0.00000466	1.0000
27	6	26	6	684.12381400	684.12383688	-0.00002288	1.0000
27	7	26	7	684.12607929	684.12610113	-0.00002184	1.0000
27	8	26	8	684.12869908	684.12871382	-0.00001474	1.0000
27	9	26	9	684.13153165	684.13167498	-0.00014333	1.0000
27	10	26	10	684.13502876	684.13498466	0.00004410	1.0000
27	11	26	11	684.13867412	684.13864291	0.00003121	1.0000
27	12	26	12	684.14266272	684.14264977	0.00001295	1.0000
27	13	26	13	684.14698268	684.14700531	-0.00002263	1.0000
27	14	26	14	684.15166360	684.15170960	-0.00004600	1.0000
27	15	26	15	684.15677738	684.15676269	0.00001469	1.0000
27	16	26	16	684.16215600	684.16216467	-0.00000867	0.0100
27	17	26	17	684.16789080	684.16791561	-0.00002481	1.0000

27	18	26	18	684.17402092	684.17401560	0.00000532	1.0000
27	19	26	19	684.18043850	684.18046473	-0.00002623	1.0000
27	20	26	20	684.18727494	684.18726308	0.00001186	1.0000
27	21	26	21	684.19435590	684.19441077	-0.00005487	1.0000
27	22	26	22	684.20186552	684.20190789	-0.00004237	1.0000
27	23	26	23	684.20981825	684.20975456	0.00006369	1.0000
27	24	26	24	684.21780671	684.21795088	-0.00014417	1.0000
27	25	26	25	684.22648900	684.22649699	-0.00000799	0.0100
27	26	26	26	684.23542400	684.23539300	0.00003100	0.0100
28	0	27	0	684.48948100	684.48918314	0.00029786	0.0033
28	2	27	2	684.48948100	684.48988024	-0.00039924	0.0033
28	1	27	1	684.48948100	684.48935742	0.00012358	0.0033
28	3	27	3	684.49062300	684.49075162	-0.00012862	0.0100
28	4	27	4	684.49189600	684.49197157	-0.00007557	0.0100
28	5	27	5	684.49353445	684.49354011	-0.00000566	1.0000
28	6	27	6	684.49544265	684.49545726	-0.00001461	1.0000
28	7	27	7	684.49776600	684.49772305	0.00004295	0.0100
28	8	27	8	684.50033101	684.50033751	-0.00000650	1.0000
28	9	27	9	684.50327558	684.50330069	-0.00002511	1.0000
28	10	27	10	684.50657403	684.50661261	-0.00003858	1.0000
28	11	27	11	684.51025271	684.51027334	-0.00002063	1.0000
28	12	27	12	684.51426949	684.51428292	-0.00001343	1.0000
28	13	27	13	684.51862494	684.51864142	-0.00001648	1.0000
28	14	27	14	684.52334900	684.52334890	0.00000010	1.0000
28	15	27	15	684.52837633	684.52840542	-0.00002909	1.0000
28	16	27	16	684.53384046	684.53381106	0.00002940	1.0000
28	17	27	17	684.53952336	684.53956591	-0.00004255	1.0000
28	18	27	18	684.54565585	684.54567003	-0.00001418	1.0000
28	19	27	19	684.55214965	684.55212353	0.00002612	1.0000
28	20	27	20	684.55890534	684.55892650	-0.00002116	1.0000
28	21	27	21	684.56607659	684.56607904	-0.00000245	1.0000

28	22	27	22	684.57367785	684.57358124	0.00009661	1.0000
28	23	27	23	684.58136992	684.58143323	-0.00006331	1.0000
28	24	27	24	684.58966200	684.58963511	0.00002689	0.0100
28	25	27	25	684.59820197	684.59818701	0.00001496	1.0000
28	26	27	26	684.60705851	684.60708905	-0.00003054	1.0000
28	27	27	27	684.61616200	684.61634136	-0.00017936	0.0100
29	0	28	0	684.86071600	684.86050253	0.00021347	0.0033
29	2	28	2	684.86071600	684.86120010	-0.00048410	0.0033
29	1	28	1	684.86071600	684.86067693	0.00003907	0.0033
29	3	28	3	684.86198300	684.86207207	-0.00008907	0.0100
29	4	28	4	684.86317400	684.86329285	-0.00011885	0.0050
29	5	28	5	684.86486847	684.86486245	0.00000602	1.0000
29	6	28	6	684.86676407	684.86678089	-0.00001682	1.0000
29	7	28	7	684.86904219	684.86904821	-0.00000602	1.0000
29	8	28	8	684.87163731	684.87166444	-0.00002713	1.0000
29	9	28	9	684.87456500	684.87462962	-0.00006462	0.0100
29	10	28	10	684.87796600	684.87794378	0.00002222	0.0100
29	11	28	11	684.88155523	684.88160698	-0.00005175	1.0000
29	12	28	12	684.88561740	684.88561927	-0.00000187	1.0000
29	13	28	13	684.88994692	684.88998071	-0.00003379	1.0000
29	14	28	14	684.89470484	684.89469137	0.00001347	1.0000
29	15	28	15	684.89977281	684.89975131	0.00002150	1.0000
29	16	28	16	684.90510717	684.90516060	-0.00005343	1.0000
29	17	28	17	684.91091240	684.91091933	-0.00000693	1.0000
29	18	28	18	684.91702500	684.91702758	-0.00000258	0.0100
29	19	28	19	684.92347569	684.92348544	-0.00000975	1.0000
29	20	28	20	684.93023991	684.93029300	-0.00005309	1.0000
29	21	28	21	684.93750415	684.93745037	0.00005378	1.0000
29	22	28	22	684.94478000	684.94495764	-0.00017764	0.0100
29	23	28	23	684.95284441	684.95281492	0.00002949	1.0000
29	24	28	24	684.96103202	684.96102234	0.00000968	1.0000

29	25	28	25	684.96963612	684.96958001	0.00005611	1.0000
29	26	28	26	684.97855800	684.97848806	0.00006994	0.0100
29	27	28	27	684.98771900	684.98774661	-0.00002761	0.0100
29	28	28	28	684.99730200	684.99735580	-0.00005380	0.0100
30	0	29	0	685.23188600	685.23152442	0.00036158	0.0017
30	2	29	2	685.23188600	685.23222245	-0.00033645	0.0017
30	1	29	1	685.23188600	685.23169893	0.00018707	0.0017
30	3	29	3	685.23300700	685.23309501	-0.00008801	0.0100
30	4	29	4	685.23426000	685.23431661	-0.00005661	0.0100
30	5	29	5	685.23590328	685.23588726	0.00001602	1.0000
30	6	29	6	685.23780532	685.23780700	-0.00000168	1.0000
30	7	29	7	685.24005963	685.24007584	-0.00001621	1.0000
30	8	29	8	685.24266352	685.24269383	-0.00003031	1.0000
30	9	29	9	685.24570784	685.24566100	0.00004684	1.0000
30	10	29	10	685.24898610	685.24897739	0.00000871	1.0000
30	11	29	11	685.25269998	685.25264305	0.00005693	1.0000
30	12	29	12	685.25660182	685.25665805	-0.00005623	1.0000
30	13	29	13	685.26104285	685.26102242	0.00002043	1.0000
30	14	29	14	685.26568057	685.26573624	-0.00005567	1.0000
30	15	29	15	685.27080300	685.27079958	0.00000342	1.0000
30	16	29	16	685.27621764	685.27621251	0.00000513	1.0000
30	17	29	17	685.28195344	685.28197511	-0.00002167	1.0000
30	18	29	18	685.28810771	685.28808747	0.00002024	1.0000
30	19	29	19	685.29454124	685.29454967	-0.00000843	1.0000
30	20	29	20	685.30143754	685.30136180	0.00007574	1.0000
30	21	29	21	685.30851034	685.30852398	-0.00001364	1.0000
30	22	29	22	685.31604990	685.31603629	0.00001361	1.0000
30	23	29	23	685.32390593	685.32389886	0.00000707	1.0000
30	24	29	24	685.33212665	685.33211179	0.00001486	1.0000
30	25	29	25	685.34068000	685.34067521	0.00000479	0.0100
30	26	29	26	685.34960040	685.34958923	0.00001117	1.0000

30	27	29	27	685.35884772	685.35885400	-0.00000628	1.0000
30	28	29	28	685.36846100	685.36846965	-0.00000865	0.0100
30	29	29	29	685.37841800	685.37843631	-0.00001831	0.0100
31	0	30	0	685.60253800	685.60224805	0.00028995	0.0033
31	2	30	2	685.60253800	685.60294656	-0.00040856	0.0033
31	1	30	1	685.60253800	685.60242268	0.00011532	0.0033
31	3	30	3	685.60355500	685.60381970	-0.00026470	0.0100
31	4	30	4	685.60495100	685.60504211	-0.00009111	0.0100
31	5	30	5	685.60661413	685.60661382	0.00000031	1.0000
31	6	30	6	685.60852918	685.60853484	-0.00000566	1.0000
31	7	30	7	685.61080191	685.61080520	-0.00000329	1.0000
31	8	30	8	685.61360600	685.61342494	0.00018106	0.0100
31	9	30	9	685.61633146	685.61639410	-0.00006264	1.0000
31	10	30	10	685.61966100	685.61971271	-0.00005171	1.0000
31	11	30	11	685.62335382	685.62338083	-0.00002701	1.0000
31	12	30	12	685.62735712	685.62739850	-0.00004138	1.0000
31	13	30	13	685.63173810	685.63176580	-0.00002770	1.0000
31	14	30	14	685.63649383	685.63648278	0.00001105	1.0000
31	15	30	15	685.64152368	685.64154950	-0.00002582	1.0000
31	16	30	16	685.64694547	685.64696606	-0.00002059	1.0000
31	17	30	17	685.65277056	685.65273251	0.00003805	1.0000
31	18	30	18	685.65878570	685.65884896	-0.00006326	1.0000
31	19	30	19	685.66538049	685.66531548	0.00006501	1.0000
31	20	30	20	685.67211609	685.67213217	-0.00001608	1.0000
31	21	30	21	685.67930278	685.67929913	0.00000365	1.0000
31	22	30	22	685.68681136	685.68681647	-0.00000511	1.0000
31	23	30	23	685.69469539	685.69468430	0.00001109	1.0000
31	24	30	24	685.70290394	685.70290272	0.00000122	1.0000
31	25	30	25	685.71144448	685.71147186	-0.00002738	1.0000
31	26	30	26	685.72051144	685.72039185	0.00011959	1.0000
31	27	30	27	685.72968624	685.72966281	0.00002343	1.0000

31	28	30	28	685.73930334	685.73928488	0.00001846	1.0000
31	29	30	29	685.74937200	685.74925820	0.00011380	0.0100
31	30	30	30	685.75989200	685.75958291	0.00030909	0.0100
32	0	31	0	685.97304900	685.97267275	0.00037625	0.0033
32	2	31	2	685.97304900	685.97337172	-0.00032272	0.0033
32	1	31	1	685.97304900	685.97284749	0.00020151	0.0033
32	3	31	3	685.97405800	685.97424544	-0.00018744	0.0050
32	4	31	4	685.97541500	685.97546867	-0.00005367	0.0100
32	5	31	5	685.97704106	685.97704142	-0.00000036	1.0000
32	6	31	6	685.97897253	685.97896372	0.00000881	1.0000
32	7	31	7	685.98122534	685.98123560	-0.00001026	1.0000
32	8	31	8	685.98389983	685.98385709	0.00004274	1.0000
32	9	31	9	685.98682955	685.98682822	0.00000133	1.0000
32	10	31	10	685.99017607	685.99014904	0.00002703	1.0000
32	11	31	11	685.99377717	685.99381960	-0.00004243	1.0000
32	12	31	12	685.99788524	685.99783996	0.00004528	1.0000
32	13	31	13	686.00224352	686.00221016	0.00003336	1.0000
32	14	31	14	686.00689521	686.00693028	-0.00003507	1.0000
32	15	31	15	686.01198625	686.01200038	-0.00001413	1.0000
32	16	31	16	686.01745153	686.01742054	0.00003099	1.0000
32	17	31	17	686.02317637	686.02319084	-0.00001447	1.0000
32	18	31	18	686.02938322	686.02931136	0.00007186	1.0000
32	19	31	19	686.03576794	686.03578218	-0.00001424	1.0000
32	20	31	20	686.04261210	686.04260341	0.00000869	1.0000
32	21	31	21	686.04977516	686.04977515	0.00000001	1.0000
32	22	31	22	686.05730327	686.05729749	0.00000578	1.0000
32	23	31	23	686.06516469	686.06517055	-0.00000586	1.0000
32	24	31	24	686.07337607	686.07339444	-0.00001837	1.0000
32	25	31	25	686.08199092	686.08196928	0.00002164	1.0000
32	26	31	26	686.09084711	686.09089520	-0.00004809	1.0000
32	27	31	27	686.10016690	686.10017233	-0.00000543	1.0000

32	28	31	28	686.10982099	686.10980080	0.00002019	1.0000
32	29	31	29	686.11982966	686.11978075	0.00004891	1.0000
32	30	31	30	686.13017933	686.13011233	0.00006700	1.0000
32	31	31	31	686.14070949	686.14079569	-0.00008620	1.0000
33	0	32	0	686.34305500	686.34279786	0.00025714	0.0033
33	2	32	2	686.34305500	686.34349729	-0.00044229	0.0033
33	1	32	1	686.34305500	686.34297272	0.00008228	0.0033
33	3	32	3	686.34414200	686.34437160	-0.00022960	0.0100
33	4	32	4	686.34551700	686.34559564	-0.00007864	0.0100
33	5	32	5	686.34718501	686.34716943	0.00001558	1.0000
33	6	32	6	686.34909299	686.34909301	-0.00000002	1.0000
33	7	32	7	686.35136498	686.35136639	-0.00000141	1.0000
33	8	32	8	686.35385363	686.35398962	-0.00013599	1.0000
33	9	32	9	686.35688900	686.35696272	-0.00007372	0.0100
33	10	32	10	686.36031580	686.36028574	0.00003006	1.0000
33	11	32	11	686.36393551	686.36395874	-0.00002323	1.0000
33	12	32	12	686.36796766	686.36798176	-0.00001410	1.0000
33	13	32	13	686.37231624	686.37235486	-0.00003862	1.0000
33	14	32	14	686.37708622	686.37707811	0.00000811	1.0000
33	15	32	15	686.38218077	686.38215157	0.00002920	1.0000
33	16	32	16	686.38755819	686.38757532	-0.00001713	1.0000
33	17	32	17	686.39333330	686.39334944	-0.00001614	1.0000
33	18	32	18	686.39947626	686.39947402	0.00000224	1.0000
33	19	32	19	686.40593403	686.40594913	-0.00001510	1.0000
33	20	32	20	686.41278245	686.41277488	0.00000757	1.0000
33	21	32	21	686.41996361	686.41995137	0.00001224	1.0000
33	22	32	22	686.42748240	686.42747869	0.00000371	1.0000
33	23	32	23	686.43535667	686.43535697	-0.00000030	1.0000
33	24	32	24	686.44364563	686.44358630	0.00005933	1.0000
33	25	32	25	686.45215699	686.45216683	-0.00000984	1.0000
33	26	32	26	686.46110912	686.46109866	0.00001046	1.0000

33	27	32	27	686.47034836	686.47038193	-0.00003357	1.0000
33	28	32	28	686.48013000	686.48001677	0.00011323	0.0100
33	29	32	29	686.49002391	686.49000333	0.00002058	1.0000
33	30	32	30	686.50021603	686.50034175	-0.00012572	1.0000
33	31	32	31	686.51088100	686.51103217	-0.00015117	0.0100
34	0	33	0	686.71294800	686.71262280	0.00032520	0.0017
34	1	33	1	686.71294800	686.71279777	0.00015023	0.0017
34	2	33	2	686.71294800	686.71332269	-0.00037469	0.0017
34	3	33	3	686.71416100	686.71419757	-0.00003657	0.0100
34	4	33	4	686.71537500	686.71542242	-0.00004742	0.0100
34	5	33	5	686.71702766	686.71699726	0.00003040	1.0000
34	6	33	6	686.71893055	686.71892210	0.00000845	1.0000
34	7	33	7	686.72115787	686.72119698	-0.00003911	1.0000
34	8	33	8	686.72386562	686.72382194	0.00004368	1.0000
34	9	33	9	686.72685511	686.72679700	0.00005811	1.0000
34	10	33	10	686.73008765	686.73012222	-0.00003457	1.0000
34	11	33	11	686.73382593	686.73379764	0.00002829	1.0000
34	12	33	12	686.73781232	686.73782331	-0.00001099	1.0000
34	13	33	13	686.74217537	686.74219930	-0.00002393	1.0000
34	14	33	14	686.74694023	686.74692567	0.00001456	1.0000
34	15	33	15	686.75204746	686.75200248	0.00004498	1.0000
34	16	33	16	686.75741419	686.75742981	-0.00001562	1.0000
34	17	33	17	686.76323972	686.76320774	0.00003198	1.0000
34	18	33	18	686.76926933	686.76933635	-0.00006702	1.0000
34	19	33	19	686.77581327	686.77581573	-0.00000246	1.0000
34	20	33	20	686.78264484	686.78264599	-0.00000115	1.0000
34	21	33	21	686.78982357	686.78982720	-0.00000363	1.0000
34	22	33	22	686.79735764	686.79735949	-0.00000185	1.0000
34	23	33	23	686.80529493	686.80524296	0.00005197	1.0000
34	24	33	24	686.81348336	686.81347772	0.00000564	1.0000
34	25	33	25	686.82203609	686.82206390	-0.00002781	1.0000

34	26	33	26	686.83098067	686.83100162	-0.00002095	1.0000
34	27	33	27	686.84033575	686.84029101	0.00004474	1.0000
34	28	33	28	686.85001914	686.84993220	0.00008694	1.0000
34	29	33	29	686.85980743	686.85992533	-0.00011790	1.0000
34	30	33	30	686.87020391	686.87027056	-0.00006665	1.0000
34	31	33	31	686.88095800	686.88096803	-0.00001003	0.0100
34	32	33	32	686.89206120	686.89201789	0.00004331	1.0000
34	33	33	33	686.90338598	686.90342031	-0.00003433	1.0000
35	0	34	0	687.08245200	687.08214702	0.00030498	0.0033
35	1	34	1	687.08245200	687.08232211	0.00012989	0.0033
35	2	34	2	687.08245200	687.08284738	-0.00039538	0.0033
35	3	34	3	687.08357400	687.08372283	-0.00014883	0.0100
35	4	34	4	687.08497700	687.08494848	0.00002852	0.0100
35	5	34	5	687.08655233	687.08652435	0.00002798	1.0000
35	6	34	6	687.08846779	687.08845046	0.00001733	1.0000
35	7	34	7	687.09078726	687.09072684	0.00006042	1.0000
35	8	34	8	687.09338836	687.09335352	0.00003484	1.0000
35	9	34	9	687.09632175	687.09633054	-0.00000879	1.0000
35	10	34	10	687.09967056	687.09965794	0.00001262	1.0000
35	11	34	11	687.10336036	687.10333578	0.00002458	1.0000
35	12	34	12	687.10735008	687.10736409	-0.00001401	1.0000
35	13	34	13	687.11175168	687.11174295	0.00000873	1.0000
35	14	34	14	687.11649440	687.11647242	0.00002198	1.0000
35	15	34	15	687.12153213	687.12155257	-0.00002044	1.0000
35	16	34	16	687.12701930	687.12698346	0.00003584	1.0000
35	17	34	17	687.13272800	687.13276518	-0.00003718	1.0000
35	18	34	18	687.13891236	687.13889782	0.00001454	1.0000
35	19	34	19	687.14539629	687.14538146	0.00001483	1.0000
35	20	34	20	687.15222388	687.15221619	0.00000769	1.0000
35	21	34	21	687.15939766	687.15940212	-0.00000446	1.0000
35	22	34	22	687.16698559	687.16693935	0.00004624	1.0000

35	23	34	23	687.17481227	687.17482799	-0.00001572	1.0000
35	24	34	24	687.18304698	687.18306816	-0.00002118	1.0000
35	25	34	25	687.19168835	687.19165997	0.00002838	1.0000
35	26	34	26	687.20068928	687.20060355	0.00008573	1.0000
35	27	34	27	687.20990100	687.20989903	0.00000197	1.0000
35	28	34	28	687.21951419	687.21954654	-0.00003235	1.0000
35	29	34	29	687.22948409	687.22954623	-0.00006214	1.0000
35	30	34	30	687.23992399	687.23989824	0.00002575	1.0000
35	31	34	31	687.25059868	687.25060272	-0.00000404	1.0000
35	32	34	32	687.26176378	687.26165983	0.00010395	1.0000
35	33	34	33	687.27309382	687.27306972	0.00002410	1.0000
35	34	34	34	687.28502400	687.28483256	0.00019144	0.0100
36	0	35	0	687.45178600	687.45137006	0.00041594	0.0017
36	1	35	1	687.45178600	687.45154527	0.00024073	0.0017
36	2	35	2	687.45178600	687.45207088	-0.00028488	0.0017
36	3	35	3	687.45286200	687.45294690	-0.00008490	0.0100
36	4	35	4	687.45423182	687.45417336	0.00005846	1.0000
36	5	35	5	687.45577431	687.45575026	0.00002405	1.0000
36	6	35	6	687.45767689	687.45767763	-0.00000074	1.0000
36	7	35	7	687.45999595	687.45995549	0.00004046	1.0000
36	8	35	8	687.46268688	687.46258389	0.00010299	1.0000
36	9	35	9	687.46551447	687.46556285	-0.00004838	1.0000
36	10	35	10	687.46890785	687.46889243	0.00001542	1.0000
36	11	35	11	687.47257307	687.47257267	0.00000040	1.0000
36	12	35	12	687.47660686	687.47660362	0.00000324	1.0000
36	13	35	13	687.48102032	687.48098534	0.00003498	1.0000
36	14	35	14	687.48569939	687.48571790	-0.00001851	1.0000
36	15	35	15	687.49081996	687.49080136	0.00001860	1.0000
36	16	35	16	687.49619571	687.49623581	-0.00004010	1.0000
36	17	35	17	687.50203980	687.50202131	0.00001849	1.0000
36	18	35	18	687.50815496	687.50815795	-0.00000299	1.0000

36	19	35	19	687.51464088	687.51464582	-0.00000494	1.0000
36	20	35	20	687.52149229	687.52148502	0.00000727	1.0000
36	21	35	21	687.52869294	687.52867565	0.00001729	1.0000
36	22	35	22	687.53621063	687.53621780	-0.00000717	1.0000
36	23	35	23	687.54412416	687.54411160	0.00001256	1.0000
36	24	35	24	687.55233476	687.55235714	-0.00002238	1.0000
36	25	35	25	687.56100640	687.56095456	0.00005184	1.0000
36	26	35	26	687.56988462	687.56990398	-0.00001936	1.0000
36	27	35	27	687.57919399	687.57920553	-0.00001154	1.0000
36	28	35	28	687.58891415	687.58885934	0.00005481	1.0000
36	29	35	29	687.59875139	687.59886555	-0.00011416	1.0000
36	30	35	30	687.60924723	687.60922432	0.00002291	1.0000
36	31	35	31	687.61995986	687.61993578	0.00002408	1.0000
36	32	35	32	687.63101409	687.63100010	0.00001399	1.0000
36	33	35	33	687.64242901	687.64241743	0.00001158	1.0000
36	34	35	34	687.65417651	687.65418795	-0.00001144	1.0000
36	35	35	35	687.66620600	687.66631181	-0.00010581	0.0100
37	0	36	0	687.82062200	687.82029151	0.00033049	0.0033
37	2	36	2	687.82062200	687.82099278	-0.00037078	0.0033
37	1	36	1	687.82062200	687.82046683	0.00015517	0.0033
37	3	36	3	687.82178100	687.82186938	-0.00008838	0.0100
37	4	36	4	687.82305900	687.82309663	-0.00003763	0.0100
37	5	36	5	687.82468937	687.82467455	0.00001482	1.0000
37	6	36	6	687.82658267	687.82660318	-0.00002051	1.0000
37	7	36	7	687.82877643	687.82888253	-0.00010610	1.0000
37	8	36	8	687.83151757	687.83151263	0.00000494	1.0000
37	9	36	9	687.83450808	687.83449354	0.00001454	1.0000
37	10	36	10	687.83783936	687.83782528	0.00001408	1.0000
37	11	36	11	687.84146564	687.84150791	-0.00004227	1.0000
37	12	36	12	687.84557061	687.84554148	0.00002913	1.0000
37	13	36	13	687.84990495	687.84992606	-0.00002111	1.0000

37	14	36	14	687.85470030	687.85466169	0.00003861	1.0000
37	15	36	15	687.85972042	687.85974846	-0.00002804	1.0000
37	16	36	16	687.86520199	687.86518644	0.00001555	1.0000
37	17	36	17	687.87098855	687.87097570	0.00001285	1.0000
37	18	36	18	687.87709100	687.87711634	-0.00002534	1.0000
37	19	36	19	687.88360477	687.88360843	-0.00000366	1.0000
37	20	36	20	687.89051817	687.89045208	0.00006609	1.0000
37	21	36	21	687.89764865	687.89764737	0.00000128	1.0000
37	22	36	22	687.90513087	687.90519443	-0.00006356	1.0000
37	23	36	23	687.91304807	687.91309335	-0.00004528	1.0000
37	24	36	24	687.92137868	687.92134426	0.00003442	1.0000
37	25	36	25	687.92991065	687.92994727	-0.00003662	1.0000
37	26	36	26	687.93889038	687.93890250	-0.00001212	1.0000
37	27	36	27	687.94822905	687.94821009	0.00001896	1.0000
37	28	36	28	687.95788820	687.95787017	0.00001803	1.0000
37	29	36	29	687.96787903	687.96788288	-0.00000385	1.0000
37	30	36	30	687.97824841	687.97824837	0.00000004	1.0000
37	31	36	31	687.98894876	687.98896679	-0.00001803	1.0000
37	32	36	32	687.99999672	688.00003829	-0.00004157	1.0000
37	33	36	33	688.01148285	688.01146304	0.00001981	1.0000
37	34	36	34	688.02322546	688.02324119	-0.00001573	1.0000
37	35	36	35	688.03540887	688.03537292	0.00003595	1.0000
37	36	36	36	688.04742500	688.04785841	-0.00043341	0.0100
38	0	37	0	688.18926100	688.18891103	0.00034997	0.0033
38	2	37	2	688.18926100	688.18961275	-0.00035175	0.0033
38	1	37	1	688.18926100	688.18908646	0.00017454	0.0033
38	3	37	3	688.19035600	688.19048991	-0.00013391	0.0100
38	4	37	4	688.19165700	688.19171796	-0.00006096	0.0100
38	5	37	5	688.19330435	688.19329691	0.00000744	1.0000
38	6	37	6	688.19518728	688.19522678	-0.00003950	1.0000
38	7	37	7	688.19755833	688.19750761	0.00005072	1.0000

38	8	37	8	688.20011022	688.20013942	-0.00002920	1.0000
38	9	37	9	688.20313545	688.20312225	0.00001320	1.0000
38	10	37	10	688.20643230	688.20645615	-0.00002385	1.0000
38	11	37	11	688.21015445	688.21014116	0.00001329	1.0000
38	12	37	12	688.21417774	688.21417735	0.00000039	1.0000
38	13	37	13	688.21852909	688.21856476	-0.00003567	1.0000
38	14	37	14	688.22329289	688.22330346	-0.00001057	1.0000
38	15	37	15	688.22838865	688.22839353	-0.00000488	1.0000
38	16	37	16	688.23379500	688.23383502	-0.00004002	0.0100
38	17	37	17	688.23963429	688.23962803	0.00000626	1.0000
38	18	37	18	688.24575338	688.24577264	-0.00001926	1.0000
38	19	37	19	688.25226200	688.25226893	-0.00000693	0.0050
38	20	37	20	688.25907036	688.25911701	-0.00004665	1.0000
38	21	37	21	688.26628676	688.26631696	-0.00003020	1.0000
38	22	37	22	688.27388953	688.27386890	0.00002063	1.0000
38	23	37	23	688.28181100	688.28177293	0.00003807	1.0000
38	24	37	24	688.29000168	688.29002918	-0.00002750	1.0000
38	25	37	25	688.29859842	688.29863775	-0.00003933	1.0000
38	26	37	26	688.30763146	688.30759877	0.00003269	1.0000
38	27	37	27	688.31690136	688.31691238	-0.00001102	1.0000
38	28	37	28	688.32658524	688.32657870	0.00000654	1.0000
38	29	37	29	688.33655631	688.33659789	-0.00004158	1.0000
38	30	37	30	688.34698471	688.34697008	0.00001463	1.0000