

TITLE:

A Workflow for Lipid Nanoparticle (LNP) Formulation Optimization Using Designed Mixture-Process Experiments and Self-Validated Ensemble Models (SVEM)

AUTHORS AND AFFILIATIONS:

Andrew T. Karl; Adsurgo, LLC; andrew.karl@adsurgo.com

Sean Essex; Arcturus Therapeutics; sean.essex@arcturusrx.com

James Wisnowski; Adsurgo, LLC; james.wisnowski@adsurgo.com

Heath Rushing; Adsurgo, LLC; heath.rushing@adsurgo.com

SUMMARY:

LNP formulation scientists frequently need to construct new recipes for different payloads or when changing lipids or process settings. We provide a robust approach to formulation optimization that minimizes the potential for errors in design construction and avoids the need for extensive statistical knowledge during analysis.

ABSTRACT:

This work summarizes a Quality by Design (QbD) styled approach to the optimization of lipid nanoparticle (LNP) formulations with a goal of providing an accessible workflow for scientists. The restriction in these studies that the molar ratios of the ionizable, helper, and PEG lipids add up to 100% requires modified design and analysis methods to accommodate this mixture constraint. Focusing on lipid and process factors that are commonly used in LNP design optimization, we provide pragmatic suggestions for avoiding the difficulties that traditionally arise in the analysis of experiments that include mixture factors and show how the recently developed statistical framework of self-validated ensemble models (SVEM) can simultaneously simplify the analysis of results from mixture-process experiments and improve the quality of the candidate optimal formulations. These steps are illustrated with a running example. We also present graphical tools based on the fitted model that simplify the interpretation of the results and facilitate the design of follow up studies in the form of confirmation runs or augmented designs.

INTRODUCTION:

Lipid nanoparticle (LNP) formulations for vivo gene delivery systems generally involve four constituent lipids from the categories of ionizable, helper, and PEG lipids (Dolgin 2022; Hou et al. 2021; Huang et al. 2022). Whether these lipids are being studied alone or simultaneously with other non-mixture factors, experiments for these formulations require “mixture” designs because – given a candidate formulation – increasing or decreasing the amount of any one of the lipids necessarily leads to a corresponding decrease or increase in the sum of the amount of the other three lipids. The ionizable and helper lipids are typically allowed to vary over a wider range of molar ratio, say 5-60%, than PEG, which we will vary from 1-5% in our example application in this paper.

Kauffman et al. (2015) consider LNP optimization studies using fractional factorial and definitive screening designs (Jones and Nachtsheim 2011); however, these types of designs cannot

accommodate a mixture constraint and are not typically used when mixture factors are present (Cornell 2002; Jones 2016). Instead, “optimal designs” capable of incorporating a mixture constraint are traditionally used for mixture-process experiments (Myers, Montgomery, and Anderson-Cook 2016). These designs target a user-specified function of the study factors and are only optimal (in one of a number of possible senses) if this function matches the true relationship between the study factors and responses. Note that there is a distinction in the text between “optimal designs” and “optimal formulation candidates”, with the latter simply being the best formulation identified by a model. Optimal designs come with three main disadvantages for mixture-process experiments. First, if the scientist fails to anticipate an interaction of the study factors when specifying the target model, then the resulting model will be biased and can produce inferior candidate formulations. Second, optimal designs place most of the runs on the exterior boundary of the factor space. In LNP studies, this can lead to a large number of lost runs if the particles do not form correctly at any extremes of the lipid or process settings. Third, scientists often prefer to have some experimental runs on the interior of the factor space to gain a model-independent sense of the response surface.

An alternative design principle is to target an approximate uniform coverage of the (mixture-constrained) factor space with a space filling design (Lekivetz and Jones 2015). These designs sacrifice some experimental efficiency relative to optimal designs (assuming the entire factor space leads to valid formulations) but present several benefits in tradeoff that are useful in this application. The space filling design does not make any a priori assumptions about the structure of the response surface, giving it flexibility to capture unanticipated aspects. This also streamlines the design generation because it does not require making decisions about which regression terms to add or remove as the desired run size is adjusted up or down. When some design points lead to failed formulations, space filling designs make it possible to model the failure boundary over the study factors while also supporting statistical models for the study responses over successful factor combinations. Finally, the interior coverage of the factor space allows for additional graphical exploration of the response surface.

To visualize the mixture factor subspace of a mixture-process experiment, specialized triangular “ternary plots” are used. Figure 1 motivates this usage: in the cube of points where three ingredients are each allowed to range from 0 to 1, the points that satisfy a constraint that the sum of the ingredients equals 1 are highlighted in red. The mixture constraint on the three ingredients reduces the feasible factor space to a triangle. In LNP applications with four mixture ingredients, we produce six different ternary plots to represent the factor space by plotting two lipids at a time against an “Others” axis that represents the sum of the other lipids.

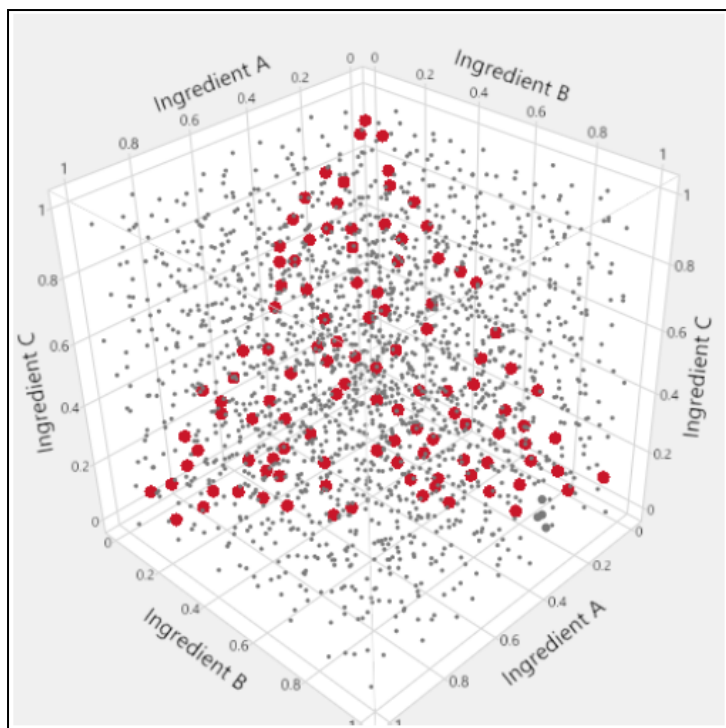


Figure 1 Triangular factor regions. In the space filling plot within the cube, the small grey dots represent formulations that are inconsistent with the mixture constraint that the three ingredients sum to 1. The larger red points lie on a triangle inscribed within the cube and represent formulations for which the mixture constraint is satisfied.

Besides the lipid mixture factors for LNP studies, there are often one or more continuous process factors such as N:P ratio, buffer concentration, or flow rate. Categorical factors may be present, such as ionizable lipid type, helper lipid type, or buffer type. The goal is to find a formulation (a mixture of lipids and settings for process factors) that maximizes some measure of potency and/or improves physiochemical characteristics such as minimizing particle size and PDI (polydispersity index), maximizing percent encapsulation, and minimizing side effects – such as bodyweight loss – in *in vivo* studies. Even when starting from a reasonable benchmark recipe, there may be interest in re-optimizing given a change in the payload or when considering changes in the process factors or lipid types.

Cornell (2002) provides a definitive text on the statistical aspects of mixture and mixture-process experiments, with Myers, Montgomery, and Anderson-Cook (2016) providing an excellent summary of the most relevant mixture design and analysis topics for optimization. However, these works can overload scientists with statistical details and with specialized terminology. We feel that modern software for the design and analysis of experiments provides a robust solution that will sufficiently support most LNP optimization problems without having to appeal to the relevant theory. While more complicated or high priority studies will still benefit from collaboration with a statistician, this workflow is meant to improve the comfort level of scientists to encourage optimization of LNP formulations without appealing to inefficient one-factor-at-a-time (OFAT, Czitrom 1999) testing or simply settling for the first formulation that satisfies specifications.

In this article, we present a workflow for utilizing statistical software for a generic LNP formulation problem, addressing design and analysis issues in the order that they will be encountered in an optimization project. Along the way, we address several common questions that arise and provide recommendations that are grounded in experience and in simulation results (Karl, Wisnowski, and Rushing 2022). Given the relatively stable geometry of the candidate factor spaces for these studies (number and ranges of mixture factors), these simulations can help identify general best practices for the design and analysis of mixture-process experiments for LNP optimization. The recently developed technique of self-validated ensemble models (SVEM, Lemkus et al. 2021) has greatly improved the otherwise fragile approach to analyzing results from mixture-process experiments, and we will show how to use this approach to use a simplified strategy that yields drastically improved formulation candidates compared to those obtained from traditional analysis methods. While the workflow is constructed in a general manner that could be followed using other software packages, JMP 17 Pro is unique in offering SVEM along with the graphical summary tools that we have found to be necessary to simplify the otherwise arcane analysis of mixture-process experiments.

SVEM employs the same linear regression model as the traditional approach, but it allows us to avoid tedious modifications that are required to fit a “full model” of candidate effects by using either a forward selection or a penalized selection (Lasso) approach, which does not need to fit the full model at once. Additionally, SVEM provides an improved “reduced model” fit that minimizes the potential for incorporating noise (process plus analytical variance) that appears in the data. It works by averaging the predicted models resulting from repeatedly reweighting the relative importance of each run in the model (Gotwalt and Ramsey 2018; Xu et al. 2020; Lemkus et al. 2021; Ramsey et al. 2021; Ramsey, Gaudard, and Levin 2021; Lemkus 2022). SVEM provides a framework for modeling mixture-process experiments that is both easier to implement than traditional single-shot regression and yields better quality optimal formulation candidates (Karl, Wisnowski, and Rushing 2022). The mathematical details of SVEM are beyond the scope of this paper and even a cursory summary beyond the relevant literature review would distract from its main advantage in this application: it allows a simple, robust, and accurate click-to-run procedure for practitioners.

The presented workflow is consistent with the Quality by Design (QbD) (Juran 1992) approach to pharmaceutical development (Yu et al. 2014). A result of the study will be an understanding of the functional relationship that link the material attributes and process parameters to critical quality attributes (CQAs, Simpson, Listak, and Hutto 2013). Daniel et al. (2022) discuss using a QbD framework specifically for RNA platform production: our workflow can be used as a tool within this framework.

In our simulated example, we want to find optimal mixtures of four constituent lipids (over three candidate ionizable lipid types), the N:P ratio, and the flow rate that maximize potency while minimizing size. The budget for the study is 23 runs (LNP batches). There is an existing control benchmark formulation that will be run within the study.

PROTOCOL:

1. Record the study purpose, responses, and factors.
 - 1.1. Write a short summary of the purpose of the experiment in a date-stamped document. Also record the choices for the design considerations in subsequent steps. For reproducibility, it is helpful to take screenshots of the design software settings.
 - 1.2. List the responses that will be measured during the experiment. Examples include potency, percent encapsulation, size, PDI, etc.
 - 1.2.1. Be careful to think of any secondary responses (e.g. downstream restrictions on physiochemical properties) that could potentially be measured simultaneously. This can prevent moving forward with a potent formulation that, for example, either causes too much bodyweight loss or will be incompatible with human injections.
 - 1.2.2. Note: when multiple assays are run for the same response on the same formulation batch, these results will be averaged prior to analysis. Unless the recipe is started from scratch, these measurements do not represent independent replicates.
 - 1.2.3. If formulation precipitation occurs or in vivo tolerability issues – such as severe body weight loss/death – are observed, then these results can be added as a column to the experimental table and modeled with binary indicators to determine the probability of this event as a function of the study factors.
 - 1.3. Determine which factors will be allowed to vary during the study and which will be held fixed.
 - 1.3.1. List the process factors that may impact one or more of the responses (Figure 2). Risk prioritization tools such as failure mode effects analyses can help select the most relevant subset for inclusion in the study (Yu et al. 2014). In most studies all lipids should be allowed to vary.

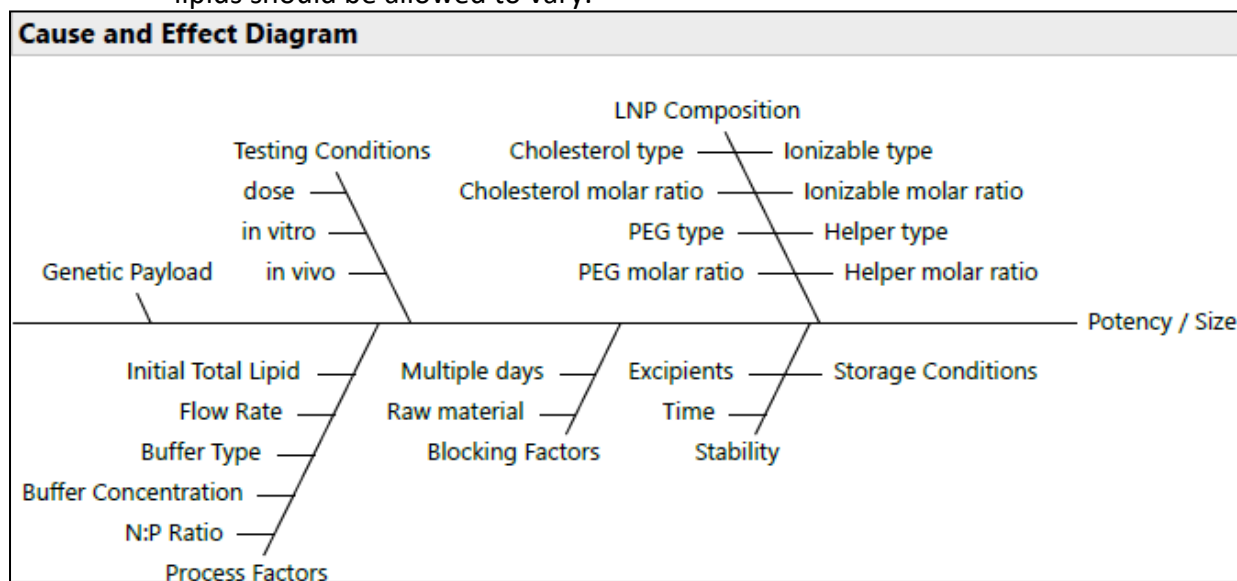


Figure 2 Cause and Effect Diagram. It is helpful to list all the factors that may affect the target responses. This figure shows common factors in an LNP formulation optimization problem.

- 1.3.2. Determine the ranges of the factors that will vary during the study.
- 1.3.3. If the study needs to be run over multiple days, then include a Day factor in the design to ensure that the factor settings are balanced across the different days. This is referred to as a “blocking” factor in experimental design literature.
- 1.3.4. Determine the relevant decimal precision for each of the factors.
 - 1.3.4.1. The factor levels listed in the design will be rounded to this level while ensuring that the mixture still adds to 100% molar ratio for the lipids.
 - 1.3.4.2. For process factors such as N:P ratio and flow rate, determine what the minimum granularity of process control exists for the factor and whether it will simplify the formulation process (and thus reduce the potential for errors) to round these factors to certain discrete values.
- 1.4. Choose the design size.
 - 1.4.1. **Minimum heuristic:** As a rule of thumb for minimum run size for continuous responses, include three runs for the optimization over each mixture factor. For four lipids, that is 12 runs. Then add two runs for each continuous process factor, and an additional run for each level of each categorical factor. For a study such as our running example with four lipid factors, two continuous and one three-way categorical process variable, this leads to a suggestion of 19 space filling runs.
 - 1.4.1.1. If it is likely that some of the runs will fail, either due to the formulations crashing out or due to measurement issues, then consider adding additional runs to compensate for the potential loss.
 - 1.4.1.2. This minimum heuristic assumes that it will be possible to perform a follow up study, if needed, in addition to performing confirmation runs for candidate optimal formulations. If it will only be possible to perform confirmation runs (e.g. for some stability studies), then it is better to budget for the number of runs obtained from the “maximum heuristic”, which is defined in the next step.
 - 1.4.2. **Maximum heuristic:** If the budget exceeds this minimum, consider including additional runs to improve the quality of the optimal formulation estimates as well as increase the power for detecting the influence of process factors. Using software for building optimal designs, the recommended run size to build a second order (main effects, two-way interactions between all effects, and quadratic effects for continuous process factors) optimal design provides a reasonable target number of runs for the space filling design employed here. For a study with four lipid factors, two continuous and one three-way categorical process variable, this leads to a suggestion of 40 space filling runs.
 - 1.4.2.1. Note: Because of diminishing marginal returns, any runs beyond this number would likely be better spent as confirmation runs or in a follow up study targeting the optimal region.
 - 1.4.2.2. Note: In JMP, the Custom Design platform can be used for this calculation.
 - 1.4.3. The design sizes derived above do not include any control benchmark runs. Add these to the total number of runs identified in the previous steps.
 - 1.4.4. Advanced note: When the primary response is binary (the measurement only

produces a yes/no response), then a larger sample size is required. It is more difficult to provide a heuristic for the required sample size for a binary response because the answer also depends on how much the response probability varies over the factor space.

2. Create the design table with a space filling design.
 - 2.1. Enter the study responses. We will illustrate by using a Potency response that is maximized and a Size response that is minimized.
 - 2.2. Optional: indicate whether each response is being maximized, minimized, or matched to a target. These settings do not affect the design and can be modified later: this is merely a convenience function.
 - 2.3. Enter the study factors and their ranges into the space filling platform of your design software (Figure 3).

| Name | Role | Values | Units |
|------------------------|-------------|----------------|-------|
| ▲ Ionizable | Mixture | 0.1 0.6 | |
| ▲ Helper | Mixture | 0.1 0.6 | |
| ▲ Cholesterol | Mixture | 0.1 0.6 | |
| ▲ PEG | Mixture | 0.01 0.05 | |
| ▼ Ionizable Lipid Type | Categorical | H101 H102 H103 | |
| ▲ N_P_ratio | Continuous | 6 14 | |
| ▲ flow rate | Continuous | 1 3 | |

Figure 3 Study factors and ranges. Screenshots of settings within experimental software are useful for reproducing the study setup.

- 2.4. Enter the number of runs for the design.
- 2.5. Optional: If using JMP, increase the average number of points per cluster from 50 to 2000 using the menu in the red triangle next to the Space Filling Design header and in the Advanced Options submenu.
- 2.6. Build the space filling design table for the chosen factors and run size. The first two runs from an example design are shown in Figure 4.

| PEG | Helper | Ionizable | Cholesterol | Ionizable Lipid Type | N_P_ratio | flow rate | Potency | Size | Notes |
|-----------|----------|-----------|-------------|----------------------|-----------|-----------|---------|------|-----------|
| 0.0419548 | 0.243548 | 0.4536548 | 0.2608424 | H103 | 8.4312589 | 2.75615 | . | . | |
| 0.01 | 0.33 | 0.33 | 0.33 | H101 | 10 | 1 | . | . | benchmark |

Figure 4 Initial output for a space filling design. Showing the first two rows of the table, settings need to be rounded to the desired precision while also making sure that the lipid amounts sum to 1. The benchmark was added to the table manually.

- 2.6.1. Add a blank “Notes” column to the table for recording information about the source of any manually created runs.
- 2.6.2. If applicable, add manual control runs to the design table. To match the target total number of runs for the experiment, either the space filling design can be built

with a smaller number of runs so that the addition of the benchmark runs will equal the target total, or the space filling design can be built with the total target number of runs and then the runs that are closes to the benchmark(s) can be replaced with the benchmark settings.

2.6.2.1. Ideally, include a replicate for the control benchmark by adding it twice to the table: this will produce a model-independent estimate of process plus analytical variance and will provide some additional graphical insight. Record the name of the benchmark in the Notes column and set the rows for the benchmark replicates to use the same color in graphs.

2.6.2.2. If any control runs lie outside of the range of the study factors, this should be recorded in the Notes column so that the row can be excluded from the subsequent analysis.

2.6.2.3. Advanced note: a “hybrid” design could be formed by creating a space filling design with most of the budgeted runs and then “augmenting” the design with the remaining runs using a D-optimal criterion (Karl, Wisnowski, and Rushing 2022). Another hybrid approach is to generate a space filling design over the mixture (lipid) and continuous (process) factors, and then afterward to add any categorical factors using an “optimal” allocation of factor levels. In JMP, this is accomplished using the “covariate” functionality of the “Custom Design” platform.

2.7. After adding any manual runs, re-shuffle the table row order. This randomization will ensure that the control runs are mixed in throughout the design rather than all appearing at the end.

2.8. Round the factor levels to the appropriate granularity and change the column formatting as needed. Format the lipid columns to appear as percentages with the desired number of decimals. After rounding the lipid ratios, ensure that they still add to 100% and make small manual changes to the factor levels if necessary. For example: the lipids may be expressed as percentages to two decimal places, the N:P ratio may be rounded to the closest multiple of 2, and flow rate may be rounded to the nearest 0.5.

2.9. Add a “Run ID” column consisting of the current date and possibly the name of the experiment concatenated with the row number of the table. The first two runs (after formatting and row randomization) now appear as Figure 5.

| Run ID | PEG | Helper | Ionizable | Cholesterol | Lipid Type | N_P_ratio | flow rate | Potency | Size | Notes |
|-----------|-------|--------|-----------|-------------|------------|-----------|-----------|---------|------|-----------|
| 02SEP22-1 | 4.19% | 24.30% | 45.36% | 26.15% | H103 | 8 | 3 | . | . | . |
| 02SEP22-2 | 1.00% | 33.00% | 33.00% | 33.00% | H101 | 10 | 1 | . | . | benchmark |

Figure 5 Formatted study table. The factor levels have been rounded and formatted and a Run ID column has been added.

2.10. Plot the design points over the lipid factors using ternary plots and examine the distribution of the runs over the process factors to ensure there are no anomalies (Figure 6 and Figure 7).

2.11. Advanced note: If any undesirable patterns appear, such as a large and uneven distribution across the factor level, the space filling design can be regenerated. If undesirable gaps appear in the ternary plots of the design, these can be addressed either with additional runs or by the hybrid approach of adding categorical factors to a

space filling design initially created over only the mixture and continuous factors.

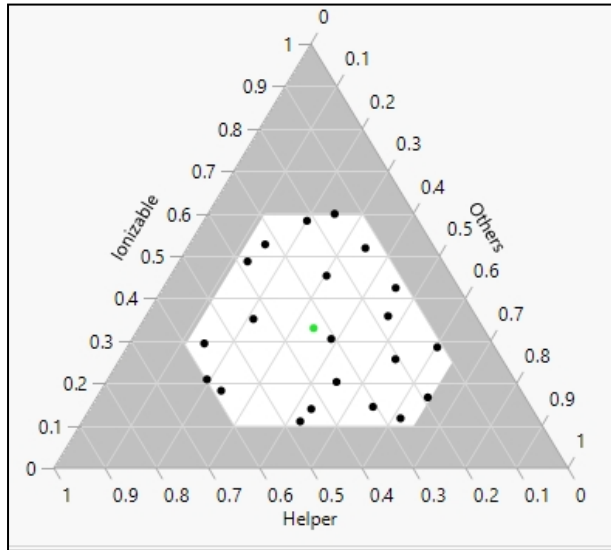


Figure 6 Design points on a ternary plot. The 23 formulations are shown as a function of their Ionizable, Helper and “Others” (Cholesterol+PEG) ratios. The green point in the center represents the benchmark.

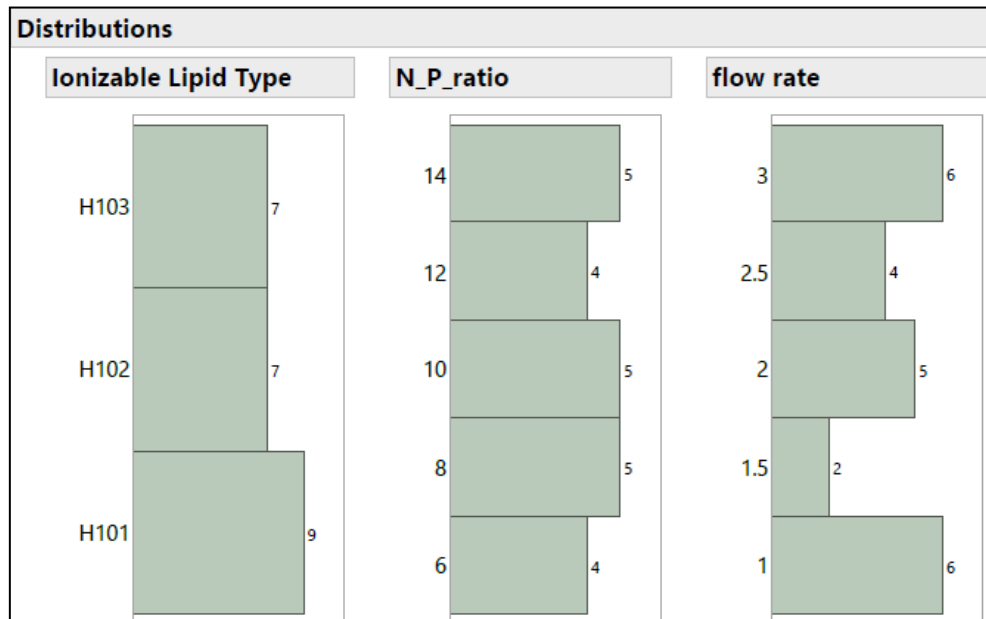


Figure 7 Distribution of non-mixture process factors in the experiment. The histograms show how the experimental runs are spaced across Ionizable Lipid Type, N:P ratio, and Flow Rate.

3. Run the experiment in the order provided by the design table.
4. Analyze the experimental results.
 - 4.1. Plot the readings and examine the distributions of the responses.
 - 4.1.1. If a replicate run was included, the relative distance between the color-coded replicate runs. Figure 8 gives an indication of the total (process and analytic) variation at the benchmark compared to the variability due to changes in the factor

settings across the factor space. Check for any concerning patterns in the plots.

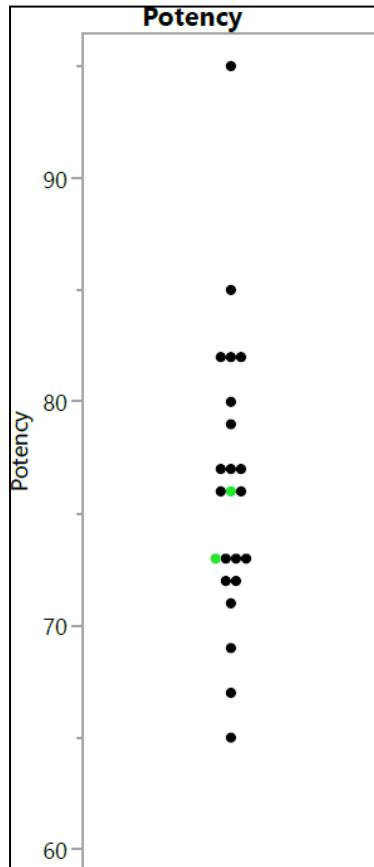


Figure 8 Observed potency readings from the experiment. The points show the potency values that were observed; the replicated benchmark runs are shown in green.

- 4.1.2. If extreme asymmetry is noted in the response distribution, then either move ahead with the normal-based regression model due to its robustness to departures from normality or consider changing the model or transforming the response.
 - 4.1.2.1. In potency or other responses with a lower bound at 0, the distribution may be so right-skewed that a few points appear to be outliers. It is possible to adjust the model either by selecting a lognormal response in a generalized regression or simply by taking a log transformation of the response prior to analysis.
 - 4.1.2.2. Percentage responses such as encapsulation sometimes show clusters of responses near the upper boundary of 100% with some left skewed readings. These can be accommodated with a beta distribution in the modeling process.
 - 4.1.2.3. Plot the runs on a ternary plot. Color the points according the responses to get a model-independent view of the behavior over the mixture factors.
- 4.2. Build an independent model for each response as a function of the study factors. Usually, the same set of candidate effects is used for each response, but a model selection technique is run independently on each response, leading to a different reduced model for each response.

- 4.2.1. Advanced note: If this is a stability study where the responses are recorded at multiple timepoints, a separate model is built for each response at each timepoint. It is often sufficient to focus on the model for the last timepoint or instead to model the differences between the last timepoint and the initial timepoint for each formulation. A joint, random-effects longitudinal modeling approach is substantially more complicated to perform and likely will not yield improved results.
- 4.2.2. Specify a “full” model of candidate effects: build a model that includes main effects for each factor along with two- and three-way interactions, quadratic and partial cubic terms in the process factors, and Scheffé cubic terms for the mixture factors (Scheffé 1958; Brown, Donev, and Bissett 2015). **Figure P8** shows some of these candidate effects. The following steps perform this action in the Fit Model platform of JMP:
- 4.2.2.1. Change the entry for **Degree** to 3 (from the default of 2). Select all of the study factors and then select **Factorial to Degree**. This adds main effects along with two- and three-way interactions.
- 4.2.2.2. Select only the mixture factors and then select **Macros > Scheffe Cubic**.
- 4.2.2.3. Select only the non-mixture factors and then select **Macros > Partial Cubic**. This adds the quadratic effects for the continuous process factors as well as the interaction of these quadratic effects with the other non-mixture factors.

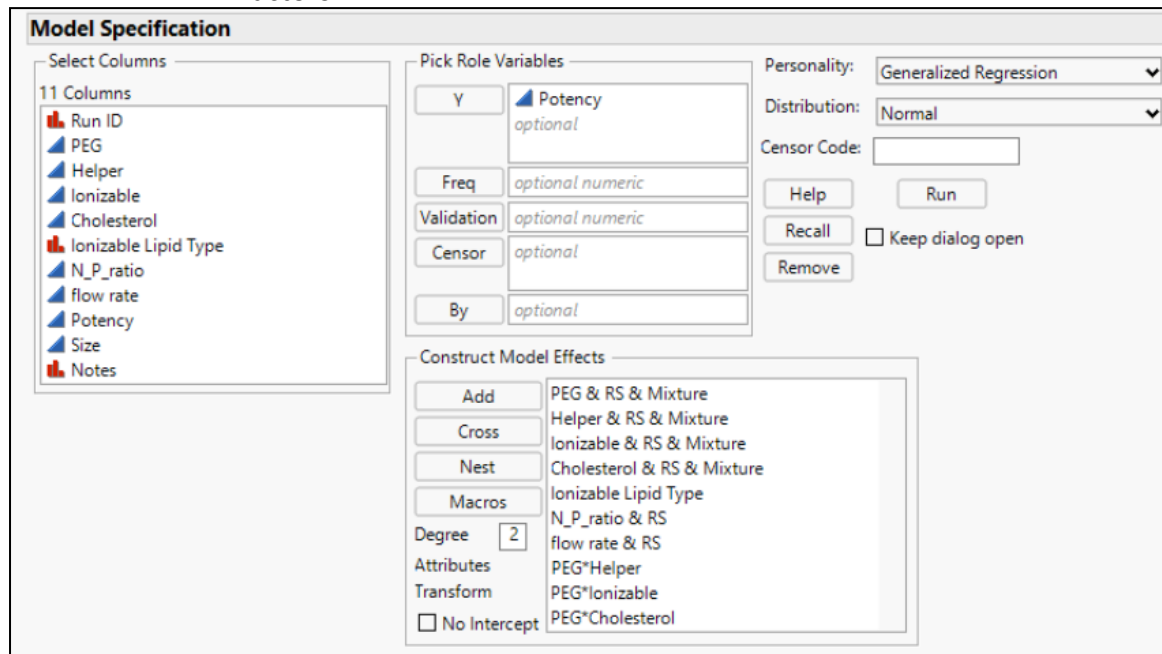


Figure 9 Software dialog for initiating the analysis. The candidate effects have been entered along with the target responses, and the No Intercept option has been unchecked.

- 4.2.3. Disable the default “No Intercept” option that analysis software typically enforces when fitting mixture effects (Figure 9). This step (to allow an intercept) is essential if employing Lasso methods. This default is in traditionally in place

because it is not possible to simultaneously fit an intercept with all of the mixture main effects without modifications (such as the SVEM approach) to the least squares regression procedure (Karl, Wisnowski, and Rushing 2022).

4.2.4. Fit a reduced model and save the prediction formula column to the data table.

4.2.4.1. Select “SVEM Forward Selection” with the appropriate response distribution, which is typically Normal, lognormal, or beta in LNP optimization problems. If unsure, then it is safest to leave the default Normal response distribution. Disable the default behavior of forcing the mixture main effects into the model. Figure 10 shows this dialog with the default behavior forcing the main effects.

Generalized Regression for Potency

Model Launch

Singularity Details

Response Distribution: Normal

Estimation Method: SVEM Forward Selection

Samples: 200

Advanced Controls

Force Terms

Forced terms are not included in the penalty.

| Force | Term |
|-------------------------------------|---------------------------------|
| <input checked="" type="checkbox"/> | Intercept |
| <input checked="" type="checkbox"/> | (PEG-0.1)/0.69 |
| <input checked="" type="checkbox"/> | (Helper-0.1)/0.69 |
| <input checked="" type="checkbox"/> | (Ionizable-0.1)/0.69 |
| <input checked="" type="checkbox"/> | (Cholesterol-0.1)/0.69 |
| <input type="checkbox"/> | Ionizable Lipid Type[H101-H103] |
| <input type="checkbox"/> | Ionizable Lipid Type[H102-H103] |
| <input type="checkbox"/> | N_P_ratio(6,14) |

Figure 10 Additional dialog for specifying SVEM options. By default, the lipid main effects are forced into the model. Because an intercept is included, we recommend unchecking these boxes in order not to force the effects.

4.2.4.2. Optional: It is also reasonable to repeat this step with the selection of “SVEM Lasso” for the Estimation Method and to see if that method leads to a different optimal recipe. If so, it can be worthwhile to run both recipes as confirmation runs to see which does best in practice (Karl, Wisnowski, and Rushing 2022).

4.2.5. For each response, plot the actual readings by the predicted readings (from the reduced model) to verify a reasonable predictive ability (Figure 11).

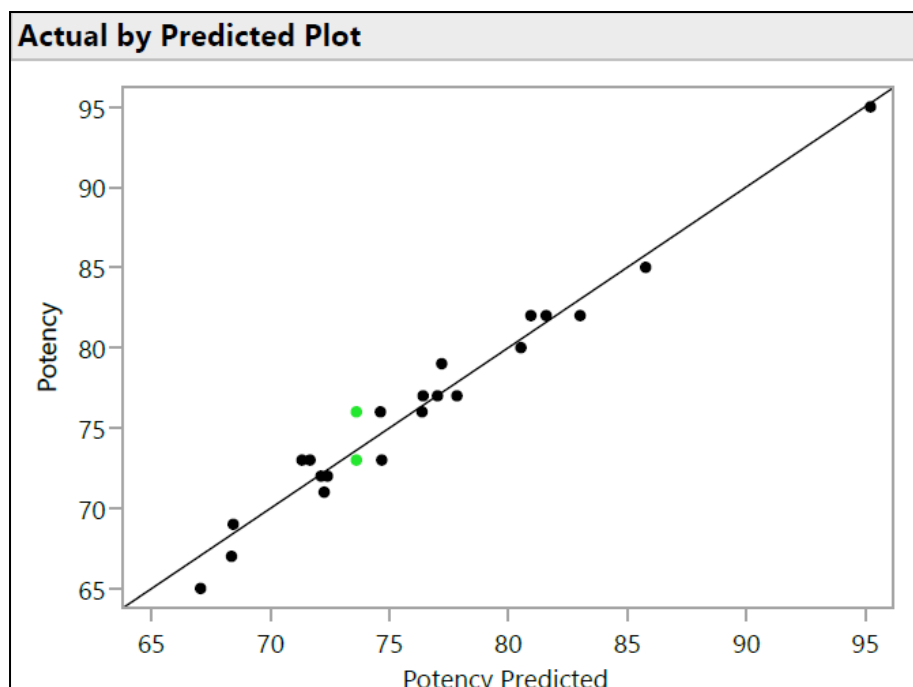


Figure 11 Actual by predicted plot. This figure plots the observed Potency against the value predicted for each formulation by the SVEM model. The correlation need not be as strong as it is in this example, but the expectation is to see at least a moderate correlation and to check for outliers.

4.3. Plot the response traces for all of the predicted response columns using a Profiler function (Figure 12).

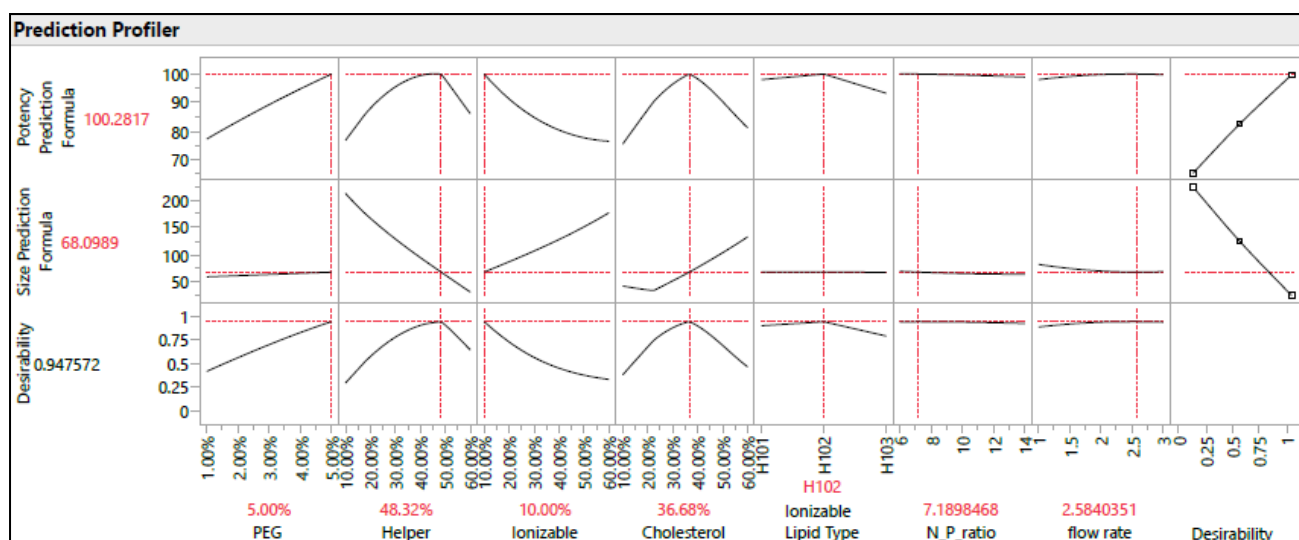


Figure 12 Prediction profiler. The top two rows of graphs show the slices of the predicted response function at the optimum formulation (as identified by the SVEM approach). The bottom row of graphs shows the weighted “desirability” of the formulation, which is a function of the last column of graphs which shows that Potency should be maximized and Size should be minimized.

4.4. Identify candidate optimal formulations.

- 4.4.1. Set the “desirability function” for each response, whether it should be maximized, minimized, or matched to a target. This also entails setting relative importance weights for each response. This might involve, for example, maximizing a primary response with importance 1.0 and minimizing a secondary response with importance 0.2. These importance weights are arbitrary relative weights between the responses; it is worthwhile to check for the sensitivity of the combined optimum to changes in these weights over a reasonable range (e.g. from equal weighting to a 1:5 weighting).
- 4.4.2. Instruct the Profiler to find the optimal factor settings that maximize the desirability function (Figure 12). NOTE: the predicted values of the responses at the optimal candidates are often overly optimistic for potency; however, the behavior at these candidate formulations will be observed directly with confirmation runs. The focus of this experimental procedure is on identifying the *location* of the optimal formulation, or the *settings* of the optimal recipe.
- 4.4.3. Record the optimal factor settings along with a note about the importance weightings used for each response.
- 4.4.4. Optional: repeat the optimization process after modifying the importance weights of the response, perhaps only optimizing the primary response or setting some of the secondary responses to have more or less importance weight (Figure 13).

| Remembered Settings | | | | | | | | | | | |
|---|------|-----------|-----------|-------------|----------------------|-----------|-----------|----------------------------|-------------------------|--------------|--|
| Setting | PEG | Helper | Ionizable | Cholesterol | Ionizable Lipid Type | N_P_ratio | flow rate | Potency Prediction Formula | Size Prediction Formula | Desirability | |
| <input type="radio"/> max potency (1.0), min size (0.2) | 0.05 | 0.4807202 | 0.1 | 0.3692798 | H102 | 7.8364699 | 2.5705497 | 100.2854 | 68.167314 | 0.989472 | |
| <input type="radio"/> max potency (1.0), min size (0.2), force H103 | 0.05 | 0.4501467 | 0.1 | 0.3998533 | H103 | 6.0139362 | 2.5562654 | 94.227002 | 78.808309 | 0.848514 | |
| <input type="radio"/> max potency | 0.05 | 0.4187891 | 0.1 | 0.4312109 | H102 | 13.89256 | 1 | 101.47714 | 100.59511 | 0.995353 | |

Figure 13 Three optimal formulation candidates from SVEM-Forward Selection. Changing the relative importance weighting of the responses can lead to different optimal formulations.

- 4.4.5. Optional: for categorical factors such as ionizable lipid type, find the conditionally optimal formulations for each factor level.
 - 4.4.6. Optional: if there is sufficient budget for an extra confirmation run, redo the analysis for each response with one more analysis technique to generate an additional joint-optimum candidate. It may be worthwhile to use the same “SVEM Forward Selection” procedure, but this time retain the default settings to build the model with “No Intercept” and to force the main effects into the model. NOTE: SVEM Lasso should never be run with the default “No Intercept” setting (Karl, Wisnowski, and Rushing 2022).
- #### 4.5. Produce graphical summaries of the optimal regions of the factor space.
- 4.5.1. Create a data table that contains 50,000 rows that are populated with randomly generated factors settings within the allowed factor space along with the predicted values from the reduced models for each of the responses and the joint desirability function. In the JMP Profiler, this is the “Output Random Table” option, and the resulting “desirability” column depends on the importance weights for the

responses that are in place when the “Output Random Table” option is selected.

4.5.2. Repeatedly produce the following graphics, changing the color-coding scheme after each iteration to produce graphs that are color coded according to the raw predicted responses of each of the responses as well as the percentile of the joint desirability function):

4.5.2.1. Ternary plots for the four lipid factors with all 50k points plotted. This will show the best and worst performing regions with respect to the lipid factors. Figure 14 shows the percentiles of the joint desirability when considering maximizing Potency (importance=1) and minimizing Size (importance=0.2), while averaging over any factors that are not shown on the ternary plot axes. Figure 15 shows the raw predicted size. NOTE: it is also reasonable to break down these graphs conditionally on other factors, such as creating a different set of ternary plots for each ionizable lipid type.

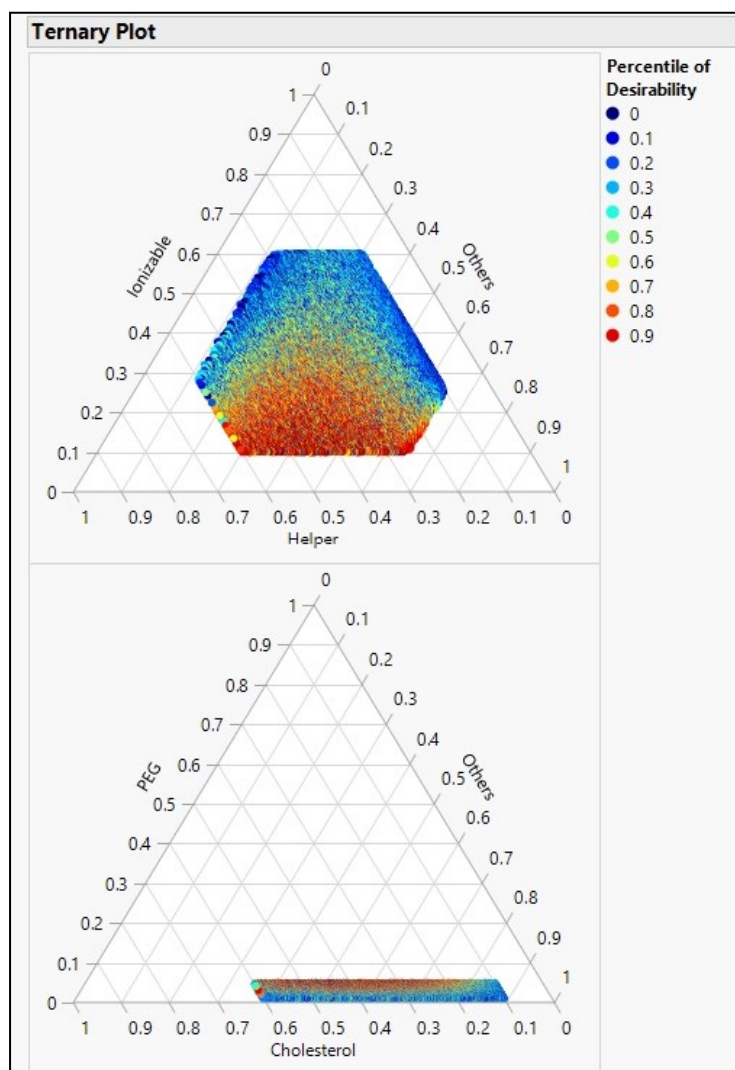


Figure 14 Ternary plots for the percentile of desirability. Where the desirability is set with importance weight of 1.0 for maximizing Potency and 0.2 for minimizing size, these plots show that the optimal region of formulations consists of lower percentages of ionizable lipid and higher percentages of PEG.

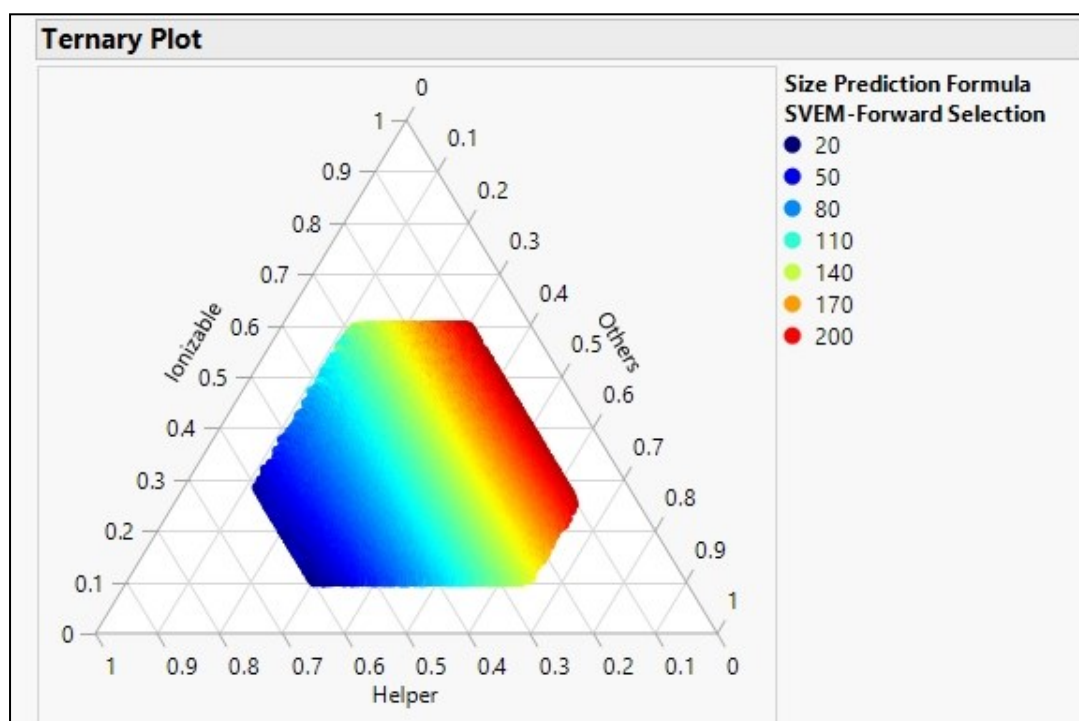


Figure 15 Ternary plot for the predicted Size. Size is minimized with higher percentages of helper lipid and maximized with lower percentages of helper.

4.5.2.2. Similarly, plot the 50k points against the non-mixture process factors, either singly or jointly. For example, plotting the 50k color-coded points with the horizontal axis set to "N:P ratio" and the vertical axis set to "flow rate" will show the jointly optimal region with respect to these factors (while averaging over all of the other factor settings). NOTE: it is also reasonable to break down these graphs conditionally on other factors, such as producing a different N:P by flow rate plot for each ionizable lipid type. Figure 16 shows the joint desirability of all of the formulations that can be formed with each of the three ionizable lipid types. The most desirable formulations use H102, with H101 providing some potentially competitive alternatives.

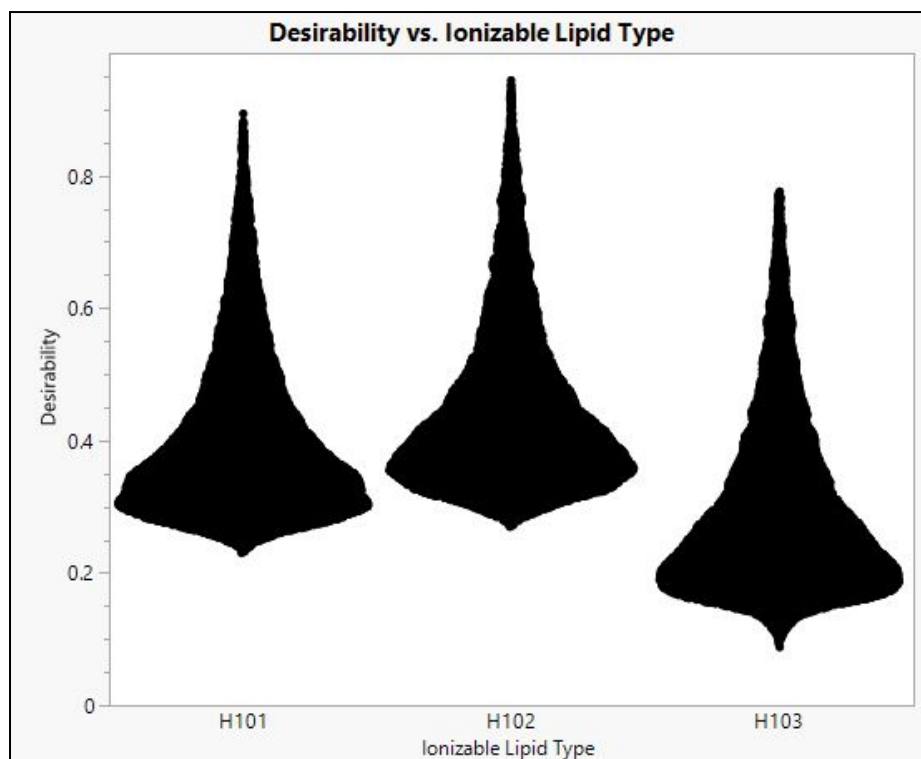


Figure 16 Violin plots for the desirability of formulations involving the three different ionizable lipid types. Each point represents a unique formulation from throughout the allowed factor space. The peaks of these distributions are the maximal values that are calculated analytically with the prediction profiler. H102 has the largest peak and thus produces the optimal formulation.

5. Plan and perform the confirmation runs.

5.1. Construct a table that lists the previously identified optimal candidates (Figure 17).

Note that the True Potency and True Size values in this table are filled in using the simulated generating functions: in reality, these will be approximated by formulating and then measuring the performance of these recipes.

| Optimal Candidate | Generating Model | PEG | Helper | Ionizable | Cholesterol | Ionizable Lipid Type | N.P.ratio | flow rate | True Potency | True Size |
|--|------------------------|-------|--------|-----------|-------------|----------------------|-----------|-----------|--------------|-----------|
| max potency (weight=1.0), min size (weight=0.2) | SVEM-LASSO | 5.00% | 47.50% | 10.00% | 37.50% | H102 | 15 | 1 | 97.4 | 71.6 |
| max potency (size ignored) | SVEM-LASSO | 5.00% | 42.35% | 10.00% | 42.65% | H102 | 15 | 1 | 98.2 | 83.8 |
| max potency (weight=1.0), min size (weight=1.0) | SVEM-LASSO | 5.00% | 51.92% | 10.00% | 33.08% | H102 | 14.2 | 3 | 95.7 | 62.6 |
| max potency (weight=1.0), min size (weight=0.2) | SVEM-Forward Selection | 5.00% | 49.25% | 10.00% | 35.75% | H102 | 14.6 | 3 | 96.8 | 67.9 |
| max potency (size ignored) | SVEM-Forward Selection | 5.00% | 42.59% | 10.00% | 42.41% | H102 | 14.9 | 1.1 | 98.2 | 83.2 |
| max potency (weight=1.0), min size (weight=1.0) | SVEM-Forward Selection | 5.00% | 52.47% | 10.00% | 32.53% | H102 | 15 | 3 | 95.4 | 61.5 |
| benchmark control | | 1.00% | 33.00% | 33.00% | 33.00% | H101 | 10 | 1 | 75.1 | 111.5 |
| H101 best: max potency (weight=1.0), min size (weight=0.2) | SVEM-Forward Selection | 5.00% | 44.33% | 10.00% | 40.67% | H101 | 6.2 | 3 | 96.9 | 78.9 |
| H103 best: max potency (weight=1.0), min size (weight=0.2) | SVEM-Forward Selection | 5.00% | 45.53% | 10.00% | 39.47% | H103 | 14.9 | 1.2 | 92.7 | 76.1 |
| best run from first experiment | | 4.74% | 42.97% | 13.96% | 38.33% | H101 | 6 | 2.5 | 93.7 | 82.2 |

Figure 17 Table of ten optimal candidates to be run as confirmation runs. The True Potency and True Size have been filled in from the simulation generating functions (without any simulated process or analytical variation).

5.1.1. If possible, include the benchmark control with the set of candidate runs.

5.1.2. If any of the formulations from the experiment proved to yield desirable results, perhaps by outperforming the benchmark, then choose the best of these formulations to carry forward and to re-run at this phase along with the optimal candidate(s) identified by the model.

- 5.1.3. Optional: Depending on the priority of the study, for example if a benchmark is going to be replaced with one of the new candidates, consider running replicates for each confirmation run, meaning that each formulation is created from scratch twice and analyzed. The average of the two results from each candidate can be used for ranking, also noting if any of the candidates seem to have a particularly large range between response values across the two replicates.
 - 5.1.4. If required, down-select from the identified candidates to match the experimental budget or to eliminate redundant candidates.
- 5.2. Perform the confirmation runs. The next step of the workflow provides instructions for constructing a follow-up study along with these runs, if desired.
- 5.3. Compare the performance of the candidate optimal formulations.
- 6. Optional: Design a follow-up study to be run concurrently with the confirmation runs
 - 6.1. Consider if one of the following is true. If not, simply perform the confirmation runs.
 - 6.1.1. The optimal formulation may lie along one of the factor boundaries, and a second experiment is desired to expand at least one of the factor ranges.
 - 6.1.2. The initial experiment used a relatively small run size and there is desire to “zoom in” on the identified optimal region with additional runs and an updated analysis.
 - 6.1.3. An additional factor is being added. This may be a level of a categorical factor (e.g. an additional ionizable lipid) or a factor that was held fixed in the initial study (e.g. buffer concentration).
 - 6.2. Build the additional experimental runs.
 - 6.2.1. Specify the new factor region in the space filling design software. This region should at least partially overlap with the region included in the initial study. If it is disconnected from the initial region, then a new study should be designed from scratch.
 - 6.2.2. Create the new experimental runs with a space filling design. Advanced users may prefer to use an optimal design instead, depending on the situation.
 - 6.2.3. After the additional space filling runs are generated, manually add two or three runs from the original experiment that are within the new factor space. These will be used to estimate any shift in the response means between blocks. Place these runs at random locations within the run-order of the experiment.
 - 6.2.4. Combine the confirmation runs and the new space filling runs into a single table and then randomize the run order.
 - 6.2.5. Advanced note: if there are a large number of confirmation runs, it is possible to account for the coverage they provide to the new factor space with a hybrid design approach that combines the confirmation runs with a number of space filling runs and then completing the design by augmenting this collection with a small number of runs from an optimal design approach. This is not a necessary step, and while it can lead to some gains in efficiency, it requires several extra steps that are beyond the scope of this workflow.
 - 6.3. Formulate the new recipes and collect the results.
 - 6.4. Concatenate the new experimental runs and results to the data table that contains the runs from the original experiment. Add a column to the data table indicating which

- experiment generated each result.
- 6.5. Analyze the results of the new experiment.
 - 6.5.1. Include the experiment-ID column as a term in the model as a blocking factor. Do not allow this term to interact with the study factors.
 - 6.5.2. Analyze the results of the new experimental results jointly with the results from the initial study. Ensure that the column properties (metadata) for each factor specifies the combined range for the factor over both studies. Follow the earlier instructions to produce updated candidate optimal formulations as well as graphical summaries.
 - 6.5.3. As a method of validation, also analyze the results from the new experiment alone, without the results from the initial experiment.
 - 6.5.3.1. Check that the optimal formulations identified by these models are in the proximity of those identified by the joint analysis.
 - 6.5.3.2. Check the graphical summaries to ensure that both the joint analysis and the individual analysis of the new experimental results report similar behavior in the response surfaces.
 - 6.5.3.3. Compare the findings of both the joint and the individual analysis of the new results to the original analysis of the initial experimental results for consistency.
 - 6.5.3.4. If any practically significant contradictions in conclusions appear then check for data entry or experimental errors. If no errors are discovered, consider any scientific explanations for the observed behavior and possibly consider an additional round of experimentation to resolve the ambiguity.
 7. Record the final scientific conclusions of the study.
 - 7.1. If the benchmark control will change to a newly identified recipe, then record the new setting and make a note about which design and analysis files record its genesis.
 - 7.2. Retain all experimental tables and analysis summaries, ideally with date-stamped file names, for future reference.

REPRESENTATIVE RESULTS:

Within the framework of the example experiment used throughout the protocol, we can repeat the simulation many times for different run sizes. A JMP script is provided for this purpose in the supplementary material. Specifically, we generate a space filling design and populate the response columns with values from our generator functions, plus noise representing analytical and process variation. We fit these simulated responses with each of the analysis techniques to produce a corresponding candidate optimal recipe. We compare the candidates from each analysis method to the value of the true optimum from the generating functions.

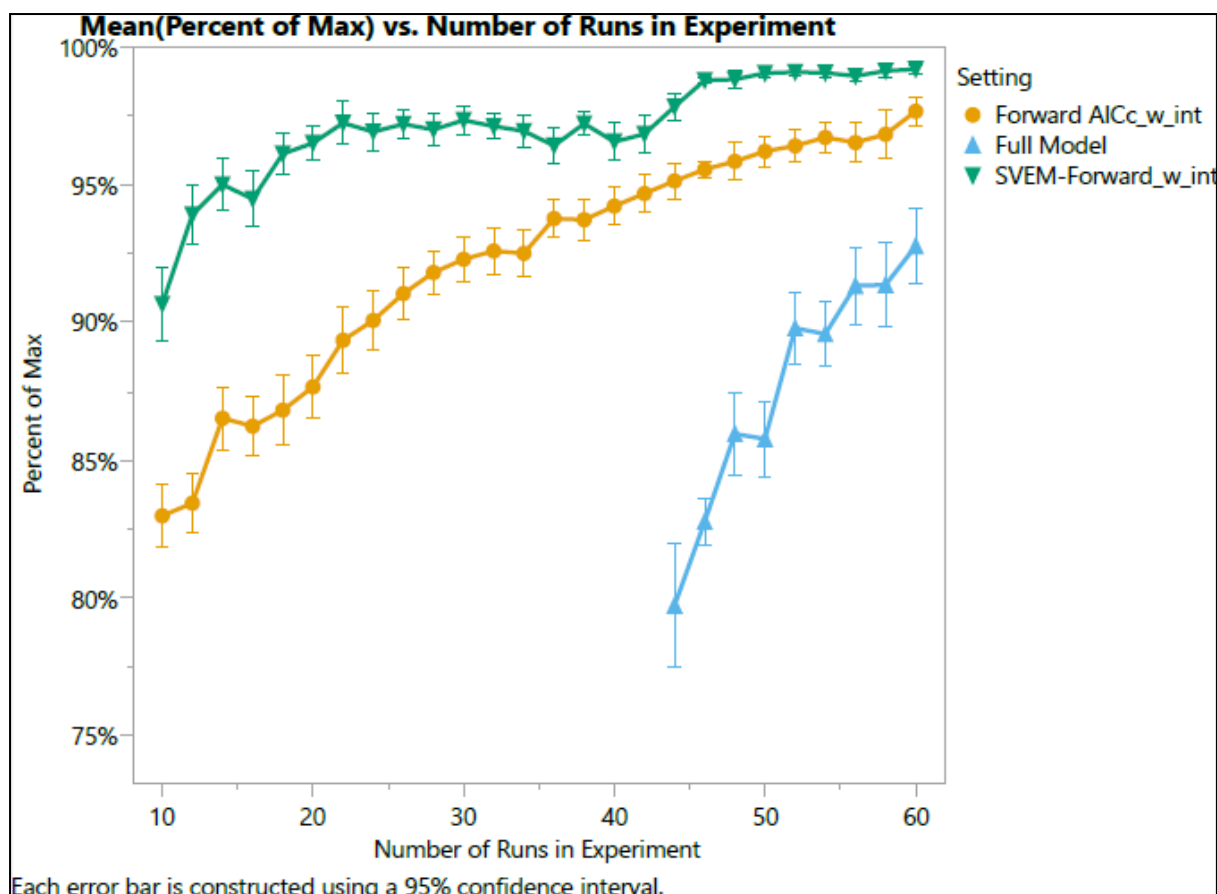


Figure 18 Quality of optimal formulation as a function of experimental size and statistical model. The vertical axis represents the percentage of theoretical maximum desirability and the horizontal axis represents the size of the space filling design. The blue line (triangles) represents the full model (without any elimination of statistically insignificant effects), the amber (circles) line represents the traditional AICc-based forward selection model (with an intercept and without forcing mixture main effects), and the green line (upside down triangles) represents the SVEM-based forward selection model.

Figure 18 shows the average percent of the maximal response achieved by each of three analysis methods using space filling designs of size given on the horizontal axis. The full model, which includes all candidate effects and no model reduction based on the statistical significance of those effects, performs the worst. Much of the additional work that traditionally goes into fitting regression models for mixture-process experiments goes into modifications (removing the intercept, forcing the mixture main effects, precluding the use of pure quadratic mixture effects, etc) that are required to fit this full model, and from this perspective those procedures are unnecessary. Furthermore, this model cannot be fit until the design size reaches the number of effects in the model. At smaller experimental sizes, we can fit the traditional forward selection method, which always outperforms the full model with respect to performance of the optimal candidate formulation for a given experiment size. Likewise, making the SVEM modification to this forward selection approach further improves the performance of the optimal candidates. In this plot, we see that using SVEM-Forward Selection to analyze a 24-run space filling experiment leads to the same average quality that would normally require 50 runs when analyzed with a traditional forward selection (targeting minimum AICc) model.

DISCUSSION:

Modern software for the design and analysis of mixture-process experiments makes it possible for scientists to improve their lipid nanoparticle formulations in a structured workflow that avoids inefficient one factor at a time (OFAT) experimentation. The recently developed SVEM modeling approach eliminates many of the arcane regression modifications and model reduction strategies that may have previously distracted scientists with extraneous statistical considerations.

The space filling designs are easier to construct – with less potential for user error – than optimal mixture-process designs. Once the results are collected, the SVEM analysis framework offers an approach that is both easier to implement and tends to produce better formulation candidates than traditional modeling approaches. Furthermore, the graphical analyses that are based on the prediction formulas for each response are easily interpretable by scientists, giving a clear summary of the marginal behavior of the response over individual factors as well as small groups of factors without requiring interpretation of highly correlated parameter estimates from a regression model.

After working through this protocol, scientists may benefit from additional reading on how optimal designs may be used in place of space filling designs, such as described by Chapter 6 of Goos and Jones (2011). Especially for the follow-up studies where there is less concern about failures along the mixture boundaries, these designs can be more efficient than space filling designs specified in the protocol. However, they require some extra decisions during setup that may frustrate inexperienced users and discourage them from using designed experiments. Furthermore, there will be occasional studies with additional complexity that will require the help of a statistician for design and analysis. Situations where the study is extremely high priority, the run size is more limited than usual, the responses are binary, there are a large number of categorical factors or levels of a single categorical factor, or where there are additional constraints on the factor space may be approached differently by the statistician, such as using optimal or hybrid (Rushing 2020) designs or adding additional structure to the design. However, this protocol has been developed over the past few years in the process of running dozens of LNP formulation optimization experiments, and we believe it offers a robust approach that will work successfully in most cases.

SUPPLEMENTARY FILES:

Available upon request from the corresponding author.

- **13DEC2022 Summary.docx** – An example document showing a record of the study purpose, responses, factors, and number of runs.
- **23_run_simulated_experiment.jmp** – A JMP file containing the simulated experiment with results. Analysis scripts that can be used with JMP 17 Pro are attached.
- **23_run_simulated_experiment.xlsx** – An Excel file containing the simulated experiment with results. For readers who do not have access to JMP.
- **mixture simulation 20DEC22.jsl** – A JMP 17 Pro script that simulates LNP formulation

experiments and records the performance of various analysis methods, including the SVEM-Forward Selection (no intercept) approach used in this workflow.

DISCLOSURES:

None.

REFERENCES:

- Brown, L., A. N. Donev, and A. C. Bissett. 2015. "General Blending Models for Data from Mixture Experiments." *Technometrics* 57 (4): 449–56.
- Cornell, John. 2002. "Experiments with Mixtures: Designs, Models, and the Analysis of Mixture Data (Wiley Series in Probability and Statistics)."
- Czitrom, Veronica. 1999. "One-Factor-at-a-Time versus Designed Experiments." *American Statistician* 53 (2): 126–31. <https://doi.org/10.1080/00031305.1999.10474445>.
- Daniel, Simon, Zoltán Kis, Cleo Kontoravdi, and Nilay Shah. 2022. "Quality by Design for Enabling RNA Platform Production Processes." *Trends in Biotechnology* 40 (10): 1213–28. <https://doi.org/10.1016/J.TIBTECH.2022.03.012>.
- Dolgin, Elie. 2022. "Better Lipids to Power next Generation of mRNA Vaccines." *Science* 376 (6594): 680–81. <https://doi.org/10.1126/SCIENCE.ADC9449>.
- Goos, Peter, and Bradley Jones. 2011. "Optimal Design of Experiments: A Case Study Approach - Peter Goos, Bradley Jones - Google Books." John Wiley & Sons, Ltd.
- Gotwalt, Chris, and Philip Ramsey. 2018. "Model Validation Strategies for Designed Experiments Using Bootstrapping Techniques With Applications to Biopharmaceuticals." In *JMP Discovery Conference*. <https://community.jmp.com/t5/Discovery-Summit-2018/Model-Validation-Strategies-for-Designed-Experiments-Using/ta-p/73730>.
- Hou, Xucheng, Tal Zaks, Robert Langer, and Yizhou Dong. 2021. "Lipid Nanoparticles for mRNA Delivery." *Nature Reviews Materials* 2021 6:12 6 (12): 1078–94. <https://doi.org/10.1038/s41578-021-00358-0>.
- Huang, Xiangang, Na Kong, Xingcai Zhang, Yihai Cao, Robert Langer, and Wei Tao. 2022. "The Landscape of mRNA Nanomedicine." *Nature Medicine* 2022 28:11 28 (11): 2273–87. <https://doi.org/10.1038/s41591-022-02061-1>.
- Jones, Bradley. 2016. "Proper and Improper Use of Definitive Screening Designs (DSDs) - JMP User Community." 2016. <https://community.jmp.com/t5/JMP-Blog/Proper-and-improper-use-of-Definitive-Screening-Designs-DSDs/bc-p/546773>.
- Jones, Bradley, and Christopher J. Nachtsheim. 2011. "A Class of Three-Level Designs for Definitive Screening in the Presence of Second-Order Effects." <https://doi.org/10.1080/00224065.2011.11917841> 43 (1): 1–15. <https://doi.org/10.1080/00224065.2011.11917841>.
- Juran, J. M. (Joseph M.). 1992. "Juran on Quality by Design : The New Steps for Planning Quality into Goods and Services," 538.
- Karl, Andrew, James Wisnowski, and Heath Rushing. 2022. "JMP Pro 17 Remedies for Practical Struggles with Mixture Experiments." In *JMP Discovery Conference*. Cary, NC. <https://community.jmp.com/t5/Discovery-Summit-Americas-2022/JMP-Pro-17-Remedies-for-Practical-Struggles-with-Mixture/ta-p/505844>.
- Kauffman, Kevin J., J. Robert Dorkin, Jung H. Yang, Michael W. Heartlein, Frank Derosa, Faryal F. Mir, Owen S. Fenton, and Daniel G. Anderson. 2015. "Optimization of Lipid Nanoparticle Formulations for mRNA Delivery in Vivo with Fractional Factorial and Definitive Screening Designs." *Nano Letters* 15 (11): 7300–7306. https://doi.org/10.1021/ACS.NANOLETT.5B02497/SUPPL_FILE/NL5B02497_SI_001.PDF.

- Lekivetz, Ryan, and Bradley Jones. 2015. "Fast Flexible Space-Filling Designs for Nonrectangular Regions." *Quality and Reliability Engineering International* 31 (5): 829–37. <https://doi.org/10.1002/QRE.1640>.
- Lemkus, Trent. 2022. "Self-Validated Ensemble Modelling." University of New Hampshire.
- Lemkus, Trent, Christopher Gotwalt, Philip Ramsey, and Maria L. Weese. 2021. "Self-Validated Ensemble Models for Design of Experiments." *Chemometrics and Intelligent Laboratory Systems* 219: 104439. <https://arxiv.org/abs/2103.09303v4>.
- Myers, Raymond, Douglas Montgomery, and Christine Anderson-Cook. 2016. "Response Surface Methodology: Process and Product Optimization Using ... -," 894. <https://www.wiley.com/en-us/Response+Surface+Methodology%3A+Process+and+Product+Optimization+Using+Designed+Experiments%2C+4th+Edition-p-9781118916018>.
- Ramsey, Philip, Marie Gaudard, and Wayne Levin. 2021. "Accelerating Innovation with Space Filling Mixture Designs, Neural Networks and SVEM." In *JMP Discovery Conference*. <https://community.jmp.com/t5/Discovery-Summit-Americas-2021/Accelerating-Innovation-with-Space-Filling-Mixture-Designs/ta-p/398669>.
- Ramsey, Philip, Wayne Levin, Trent Lemkus, and Christopher Gotwalt. 2021. "SVEM: A Paradigm Shift in Design and Analysis of Experiments." In *JMP Discovery Conference - Europe*. <https://community.jmp.com/t5/Discovery-Summit-Europe-2021/SVEM-A-Paradigm-Shift-in-Design-and-Analysis-of-Experiments-2021/ta-p/349244>.
- Rushing, Heath. 2020. "DOE Gumbo: How Hybrid and Augmenting Designs Can Lead To More Effective Design Choices." In *JMP Discovery Conference*.
- Scheffé, Henry. 1958. "Experiments with Mixtures." *Journal of the Royal Statistical Society: Series B (Methodological)* 20 (2): 344–60. <https://doi.org/10.1111/J.2517-6161.1958.TB00299.X>.
- Simpson, James R, Charles M Listak, and Gregory T Hutto. 2013. "Guidelines for Planning and Evidence for Assessing a Well-Designed Experiment." *Quality Engineering* 25 (4): 333–55. <https://doi.org/10.1080/08982112.2013.803574>.
- Xu, Li, Chris Gotwalt, Yili Hong, Caleb B King, and William Q Meeker. 2020. "Applications of the Fractional-Random-Weight Bootstrap View Supplementary Material Applications of the Fractional-Random-Weight Bootstrap." <https://doi.org/10.1080/00031305.2020.1731599>.
- Yu, Lawrence X., Gregory Amidon, Mansoor A. Khan, Stephen W. Hoag, James Polli, G. K. Raju, and Janet Woodcock. 2014. "Understanding Pharmaceutical Quality by Design." *The AAPS Journal* 16 (4): 771. <https://doi.org/10.1208/S12248-014-9598-3>.