
RAD-PHI2: INSTRUCTION TUNING PHI-2 FOR RADIOLOGY

A PREPRINT

Mercy Ranjit*
Microsoft Research, India
meranjit@microsoft.com

Gopinath Ganapathy
Department of Computer Science
Bharathidasan University, India
gganapathy@bdu.ac.in

Shaury Srivastav
Microsoft Research, India
t-ssrivastav@microsoft.com

Tanuja Ganu
Microsoft Research, India
taganu@microsoft.com

Srujana Oruganti
Microsoft Research, India
sruvs@microsoft.com

March 18, 2024

ABSTRACT

Small Language Models (SLMs) have shown remarkable performance in general domain language understanding, reasoning and coding tasks, but their capabilities in the medical domain, particularly concerning radiology text, is less explored. In this study, we investigate the application of SLMs for general radiology knowledge specifically question answering related to understanding of symptoms, radiological appearances of findings, differential diagnosis, assessing prognosis, and suggesting treatments w.r.t diseases pertaining to different organ systems. Additionally, we explore the utility of SLMs in handling text-related tasks with respect to radiology reports within AI-driven radiology workflows. We fine-tune Phi-2, a SLM with 2.7 billion parameters using high-quality educational content from Radiopaedia [Radiopaedia.org], a collaborative online radiology resource. The resulting language model, RadPhi-2-Base, exhibits the ability to address general radiology queries across various systems (e.g., chest, cardiac). Furthermore, we investigate Phi-2 for instruction tuning, enabling it to perform specific tasks. By fine-tuning Phi-2 on both general domain tasks and radiology-specific tasks related to chest X-ray reports, we create Rad-Phi2. Our empirical results reveal that Rad-Phi2 Base and Rad-Phi2 perform comparably or even outperform larger models such as Mistral-7B-Instruct-v0.2 and GPT-4 providing concise and precise answers. In summary, our work demonstrates the feasibility and effectiveness of utilizing SLMs in radiology workflows both for knowledge related queries as well as for performing specific tasks related to radiology reports, thereby opening up new avenues for enhancing the quality and efficiency of radiology practice.

Keywords Small Language Models (SLMs) · Instruction Tuning · Radiology Reports

1 Introduction

Recent advances in natural language processing (NLP) have enabled the development of powerful language models that can generate natural language texts for various purposes and domains. However, most of the existing language models are trained on general domain texts, such as news articles, books, web pages, etc., and may not be able to handle the specific vocabulary, concepts, terminology, and style of radiology texts. Hence, it is required to have domain specific models that embody knowledge for that domain and trained for specific tasks in domain specific workflows to deal with the nuances of the field. But, the size and complexity of these language models are often prohibitive for practical applications, as they require a large amount of computational resources and memory to train and run.

*The author is also affiliated as a PhD student with Bharathidasan University.

In this paper, we explore the use of small language models (SLMs) for radiology related tasks, specifically text-related tasks relating to radiology knowledge and tasks that are useful in everyday radiology workflows pertaining to radiology reports.

SLMs are language models that have a relatively small number of parameters, compared to the state-of-the-art large language models, such as GPT-4 [OpenAI et al., 2023]. SLMs, specifically Phi-2 [Gunasekar et al., 2023] is a 2.7 billion parameter model that has shown remarkable performance in general domain language understanding, reasoning, and coding tasks, but its capabilities in the medical domain, especially on radiology content and tasks specific to a radiology workflows is not explored. Moreover, Phi-2 is not an instruction tuned model and tends to produce verbose answers as it was trained from synthetic text books and currently not suitable for direct deployments in application workflows to perform specific tasks. We hypothesize that SLMs can be effectively adapted to the radiology domain by finetuning them on high-quality radiology content to bring-in the nuanced radiology knowledge and also be adapted to perform specific tasks that are useful in everyday radiology workflows. SLMs provide the advantage of having a domain specific model for easy consumption at lower computational cost and faster inference.

The Phi models [Gunasekar et al., 2023] success and its on-par performance to large language models was attributed to the high quality aka "*textbook-quality*" training data that was used for its training. The training corpus was augmented with carefully selected data, filtered based on educational value and content quality. It demonstrated that this choice of training data enabled Phi-2 to achieve outstanding performance on complex benchmarks, matching or outperforming models up to 25x larger in size. We extend this approach of Phi-2 by using high-quality and reliable educational content from Radiopaedia [Radiopaedia.org] to finetune Phi-2 for radiology knowledge. Radiopaedia is a rapidly growing collaborative online peer-reviewed open-edit radiology resource, compiled by radiologists and other health professionals from across the globe. Training from this highly valuable and reliable source can embody the model with this knowledge making it useful to answer questions on radiology accurately. To test our hypothesis, we use Phi-2, a 2.7 billion parameter SLM, as our base model, and finetune it on question-answer pairs derived from the educational articles of Radiopaedia pertaining to different systems like chest, cardiac, gastrointestinal etc. The resulting model Rad-Phi2 Base, can answer queries on radiology pertaining to these systems accurately. We evaluate the performance of Rad-Phi2-Base and compare it with larger instruction tuned models, such as Mistral-7B-Instruct-v0.2 [Jiang et al., 2023] which was finetuned on the same radiology question answer pairs and also GPT-4. Rad-Phi2 Base performs better than or comparable to the much larger counterparts while producing more crisp and to-the-point answers. It even surpassed the Retrieval Augmented Generation (RAG) approach [Lewis et al., 2021] that directly used a knowledge base of Radiopaedia articles to answer the queries.

Apart from answering queries on radiology, there are specific text-related tasks that can be useful in AI deployment scenarios of radiology workflows. Some useful tasks relating to radiology reporting workflows are writing diagnostic impression summary from radiology findings, tagging of radiology reports with appropriate abnormality labels to flag priority cases, noise free extraction of findings and impression from existing radiology reports, generating a comparison summary of findings between the current radiology report and its prior report etc. Some of these tasks were also cited in the earlier work [Bhayana, 2024] in the context of application of language models in radiology workflows. To address these useful tasks that arise in the AI deployment scenarios of radiology workflows, we create an instruction tuning dataset for tasks pertaining to radiology reports of Chest X-rays and instruction-tune Phi-2 and obtain the Rad-Phi2 model. But, as Phi-2 was not an instruction tuned model in the first place, we first instruction tune Phi-2 for general domain tasks using the Super Natural instruction tuning dataset [Yizhong and Swaroop, 2022]. The resulting model Phi2-Instruct (SupNatInst) is then further instruction tuned for tasks specific to radiology reports. The resulting model, Rad-Phi2 can perform specific tasks that are useful in AI driven radiology reporting workflows of Chest X-rays. Our results show that Rad-Phi2 performs better or on-par with a larger model like Mistral-7B-Instruct (v0.2) finetuned on the same instruction tuning dataset.

We would also like to note that while GPT-4 has reasonable radiology knowledge as noted in this work [Liu et al., 2023a], it still falls short on nuanced tasks like impression generation in comparison to domain specific models unless given with a large number of carefully curated examples in conjunction with sophisticated prompt engineering approaches [Ma et al., 2024]. SLMs enables to develop domain-specific models for these nuanced task types easily and quickly due to the low resource requirements associated with its size.

Our work demonstrates the feasibility and effectiveness of using SLMs for radiology specific tasks, and opens up new possibilities for improving the quality and efficiency of radiology practice.

Our key Contributions include:

- We demonstrate the effectiveness of SLMs in the radiology domain by training on high quality radiology content from Radiopaedia [Radiopaedia.org]. The resulting model Rad-Phi2-Base can answer general queries on radiology like symptoms, radiological appearances of findings, differential diagnosis, prognosis,

and treatments w.r.t diseases pertaining to 15 different systems like Chest, Cardiac etc. It performs better than larger models like `Mistral-7B-Instruct-v0.2` finetuned on the same dataset. It also performed better than GPT-4 and the Retrieval Augmented Generation approach that directly used Radiopaedia content as the knowledge source while answering queries.

- We construct an instruction tuning dataset of text related tasks that are useful in radiology report writing workflows of Chest X-rays using GPT-4 processing and annotations of existing datasets on Chest X-ray reports. We share the prompts for different tasks used in creating this dataset.
- We also demonstrate Instruction Tuning on SLMs to make it usable for practical tasks involved in radiology workflows. As Phi-2 is not an instruction tuned model and tends to produce verbose answers, we first perform general domain instruction tuning using the Super Natural Instructions dataset followed by radiology specific instruction tuning on tasks related to radiology reports of Chest X-rays. We call the resulting model Rad-Phi2 which is useful in performing everyday tasks in radiology report writing workflows of Chest X-rays.
- Our ablation studies gives insights on the input formats that worked the best when instruction tuning the Phi-2 model. We note that radiology specific instruction tuning performs better when preceded by general domain instruction tuning and semantic tasks like Natural Language Inference(NLI) can help with the performance of other tasks.

2 Related Work

The work, ‘Exploring the Boundaries of GPT-4 in Radiology’ [Liu et al., 2023a] evaluated GPT-4 on a diverse range of common radiology text-based tasks. It was noted that that zero-shot prompting technique was sufficient for the ‘*learn-by-instruction*’ tasks like sentence similarity, radiology NLI, and disease progression. However, a large number of carefully curated examples with advanced prompt engineering were required for ‘*learn-by-example*’ tasks like impression summary writing to perform on par with the supervised models trained for the task. While it is beneficial to have GPT-4 like instruction-tuned models that can work well on radiology content for ‘*learn-by-instruction*’ tasks, it is still beneficial to have specialized models for ‘*learn-by-example*’ kind of tasks where the required knowledge for the model needs to be fed in as carefully curated in-context examples. It is even better if a model can do both these kinds of tasks with high accuracy at lower computational and resource costs, which is what we accomplish in our work.

The work ‘Chatbots and Large Language Models in Radiology: A Practical Primer for Clinical and Research Applications’[Bhayana, 2024] highlights how LLMs are useful in typical radiology workflows, providing key relevant information from clinical notes and prior reports, thus improving radiology interpretation and clinical decision making. The paper also noted LLMs when connected to sources like the American College of Radiology performed better than radiologists in applying this information to clinical scenarios. In regular practice, radiologists and clinicians refer to external knowledge sources like Radiopaedia [Radiopaedia.org] for information like clinical symptoms, radiological appearances of findings, differential diagnosis, prognosis, and treatment plan of a disease that aids them with clinical decision making. These knowledge sources are very reliable as they are maintained by practitioners of medical practice and are actively updated. However, this also poses a need for the AI models to be updated regularly with this knowledge. While Retrieval Augmented Generation [Lewis et al., 2021] can help LLMs augment their responses with the knowledge from the external sources, it may prove costly for high usage scenarios. This calls for a need for having Small Language Models that are easy to train, deploy, retrain, and iterate.

While instruction tuning to perform multiple tasks has been well explored on general domain tasks, it is less explored on radiology specific tasks. Radiology-Llama2[Liu et al., 2023b] and Radiology-GPT[Liu et al., 2023c] performed instruction tuning for one specific task of generating clinically useful impression summary from radiological findings. However, there are a variety of tasks that can be useful in radiology report writing workflows, like extracting a clean set of findings or impressions from a noisy radiology report, tagging the reports with abnormality labels or tubes, lines and devices labels, or comparing two reports to see the differential findings between them, querying the reports for specific questions or querying for general radiology knowledge with respect to various diseases like their radiographic appearances, symptoms, differential diagnosis or treatment plans to make a clinical decision. We consider all these useful tasks for instruction tuning.

We explore the usage of Small Language Models for Radiology-specific tasks, as these models are very easy to train and deploy due to their smaller size, while being very competitive in performance to large language models. Phi-2 paper [Gunasekar et al., 2023] introduced a 2.7 billion-parameter language model that demonstrated outstanding reasoning and language understanding capabilities, showcasing state-of-the-art performance among base language models with less than 13 billion parameters. On complex benchmarks, Phi-2 matched or outperformed models up to 25x larger in size. The authors attributed this success to the high-quality textbook-like data they trained the model on. This is very applicable in the medical field, where the information rendered by the model needs to be accurate, which mandates the

Table 1: RADIOPAEDIA QUESTION ANSWERING DATA STATISTICS

Systems	Article Counts	QA Pair Counts
Chest	1741	9695
Cardiac	771	5188
Central Nervous System	2846	16557
Urogenital	705	4171
Oncology	391	2782
Breast	357	1929
Musculoskeletal	3804	24026
Not Specified	1579	10880
Hepatobiliary	470	2956
Vascular	423	2304
Gastrointestinal	1240	7342
Obstetrics	568	3157
Interventional	149	1003
Trauma	85	483
Spine	94	490
Forensic	18	105
Total	15076	93068

Count of total articles and respective question answer pairs by system.

training dataset to be noise-free and of high quality. Our work considers this need of very high-quality data for training the models, which is practically deployable in day-to-day radiology workflows.

We also specifically note that Phi-2 is a base model that has not undergone alignment through reinforcement learning from human feedback (RLHF) [Ouyang et al., 2022], nor has it been instruction finetuned. Due to this, it is noted that it often produces irrelevant or extra text following its first answer to user prompts within a single turn. This is due to its training dataset being primarily textbooks, which results in textbook-like responses. Our work addresses these issues by instruction tuning Phi-2 first for general domain tasks to make it usable for downstream applications in Radiology.

In summary, our paper addresses the need for having special instruction tuned models for radiology which has accurate radiology knowledge and also can perform useful tasks on radiology workflows which is not fully addressed in the above stated works. In addition, we also perform these tasks using SLMs that are easy to train, deploy and iterate on without compromising on quality.

3 Method

3.1 Dataset and Tasks

General Domain Instruction Tuning We performed minimal general domain instruction tuning before we performed radiology specific finetuning using the Super Natural Instructions dataset [Yizhong and Swaroop, 2022]. This dataset is a benchmark of 1,616 diverse NLP tasks including non-english tasks and their expert-written instructions. The collection covers 76 distinct task types. We used only the english specific tasks with a total of 757 tasks and 2,410,002 training records corresponding to these tasks.

Radiology Specific Instruction Tuning We create two datasets for radiology specific finetuning, one for general radiology question answering and another for instruction tuning for tasks specific to radiology reports of Chest X-rays.

For radiology question answering task, we create question answer pairs using GPT-4 from the articles of Radiopaedia pertaining to 15 different systems. We reserve the summary articles pertaining to these systems for testing. The count of question answer pairs w.r.t different systems is detailed in Table 1. This dataset was used to train the Rad-Phi2-Base model.

The instruction tuning dataset for tasks specific to radiology reports was constructed from the radiology reports of Mimic-CXR dataset [Johnson et al., 2019] that contains radiology reports pertaining to Chest X-rays. We use GPT-4 processing and also annotations from the existing datasets derived from Mimic-CXR reports to construct this dataset. The prompts associated with the different tasks in the instruction tuning dataset is defined in Table 2. This dataset was

Table 2: PROMPT DESIGN - INSTRUCTION TUNING DATASET FOR RADIOLOGY REPORTS

Task Name	Prompt
RADIOLOGY QA	<instruct>Question <output>Answer<endoftext>
EXTRACT IMPRESSION	<instruct>Given the radiology report below: FULL REPORT Output the findings from the findings section of the radiology report. Each sentence in the output should describe an observation or a finding about the image. Do not mention any changes in observations, follow-up suggestions, doctor names, or noisy details. <output>Answer<endoftext>
EXTRACT FINDINGS	<instruct>Given the radiology report below: FULL REPORT Output the impression of the radiology report. Each sentence in the output should describe an observation or a finding about the image. Do not mention any changes in observations, follow-up suggestions, doctor names, or noisy details. <output>Answer<endoftext>
CLEANUP REPORT TEXT	<instruct>Given the text from a radiology report: RADIOLOGY REPORT TEXT Update the impressions or findings such that each sentence in the output describes an impression or observation about the image. Remove any mention of change of an observation and just state its presence. Do not include any follow-up suggestions or advice, and avoid mentioning any doctor names or other noisy details. <output>Answer<endoftext>
QA COMPREHENSION	<instruct>Answer the question using the radiology report below as context: FULL REPORT Question: <output>Answer<endoftext>
QA TEMPORAL FINDINGS	<instruct>Given the below radiology report for an image and its prior report for reference: FULL REPORT PRIOR REPORT What findings are added and what findings are removed in the current radiology report for an image, compared to its reference report from before? <output>Answer<endoftext>
QA TEMPORAL PROGRESSION	<instruct>Given the radiology report below, classify the progression of a finding as improved, no change, worsened. FULL REPORT <output>Answer<endoftext>
ABNORMALITY LABELS	<instruct>Given the below radiology report: FULL REPORT What abnormality labels can be tagged to these findings? <output>Answer<endoftext>
TUBES, LINES AND DEVICES LABELS	<instruct>Given the below radiology report: FULL REPORT Identify the tubes and lines or devices that are mentioned in the radiology report above. <output>Answer<endoftext>
IMPRESSION PREDICTION	<instruct>Given the findings from a radiology report: FINDINGS Based on the above findings from a radiology report, write an impression. <output>Answer<endoftext>

Prompt Design for the Instruction Tuning dataset that was used to train Rad-Phi2.

Table 3: RADIOLOGY REPORTS INSTRUCTION TUNING DATASET

Task	Train	Test	Validate	Total
QA TEMPORAL PROGRESSION	50000	17639	7078	74717
ABNORMALITY LABELS	221035	3403	1959	226397
QA TEMPORAL FINDINGS	100000	2939	1308	104247
TUBES AND LINES DEVICES LABELS	94915	3403	1959	100277
IMPRESSION PREDICTION	208876	2523	1647	213046
EXTRACT FINDINGS	50000	3844	2103	55947
EXTRACT IMPRESSION	50000	3283	2295	55578
CLEANUP RADIOLOGY TEXT	75000	7127	4398	86525
QA COMPREHENSION	467057	9179	3878	480114
Total	1316883	53340	26625	1396848

Rad-Phi2 Instruction Tuning Dataset Statistics: Count of task wise splits for train, test and validation sets.

used to train the Rad-Phi2 model. The number of training records in the instruction tuning dataset for various tasks is detailed in Table 3.

We construct the below tasks:

EXTRACT FINDINGS: Extract the findings from a given radiology report.

EXTRACT IMPRESSION: Extract the impression section from a given radiology report.

CLEANUP RADIOLOGY REPORT TEXT: Given a radiology report text of findings or impression which includes noisy details like mention of Dr. names, followup actions etc., remove the noisy details to produce a clean radiology text of findings. We used GPT-4 processing to generate a clean set of findings from a radiology report.

QA COMPREHENSION: Given a radiology report, answer the question based on the radiology report. The question could be about the abnormalities in the report, associated location or severity of findings etc. We use the Medical-Diff VQA dataset [Hu et al., 2023] to construct the QA comprehension task.

QA TEMPORAL FINDINGS: Given a radiology report and its prior report, extract the differential set of findings between them. We use the Medical-Diff VQA dataset [Hu et al., 2023] to construct the QA temporal comprehension task from the corresponding radiology reports.

QA TEMPORAL PROGRESSION: Given a radiology report, classify the progression of an abnormality as improved, worsened or unchanged. We use the ChestImagenome dataset [Wu et al., 2021] for the progression labels.

ABNORMALITY LABELS: Given a radiology report, label the report as positive for a list of abnormalities. We tag the reports to findings/disease labels defined in the ChestImagenome dataset [Wu et al., 2021].

TUBES, LINES AND DEVICES LABELS: Given a radiology report, label the report for the presence of tubes, line and devices. We tag the reports to tubes, lines and device labels defined in the ChestImagenome dataset [Wu et al., 2021].

IMPRESSION PREDICTION: Given the findings from a radiology report, summarize the findings as a diagnostic impression.

3.2 Training Workflow

The training workflow is defined in the Fig. 1. We started with the Phi-2 model and first instruction tuned it using the general domain tasks defined in the Super Natural Instructions dataset [Yizhong and Swaroop, 2022]. We refer the resulting version as Phi2-Instruct(SupNatInst). This was further finetuned for Radiology domain to produce two different models, Rad-Phi2-Base for answering radiology specific questions pertaining to 15 different systems defined in Table 1 and Rad-Phi2 for the tasks specific to radiology reports defined in Table 3.

We use GPT-4 for processing the articles from Radiopaedia [Radiopaedia.org] into question answer pairs for creating the Radiology QA dataset that was used to train the Rad-Phi2-Base. We also use it for processing the Radiology Reports to extract a clean set of findings from the radiology reports used in the creation of the Extract Findings/Impression and Cleanup Radiology Text tasks of the instruction tuning dataset. We also use annotations from existing datasets and apply the prompts defined in Table 2 with these annotations to create the instruction tuning dataset that was used in training Rad-Phi2.

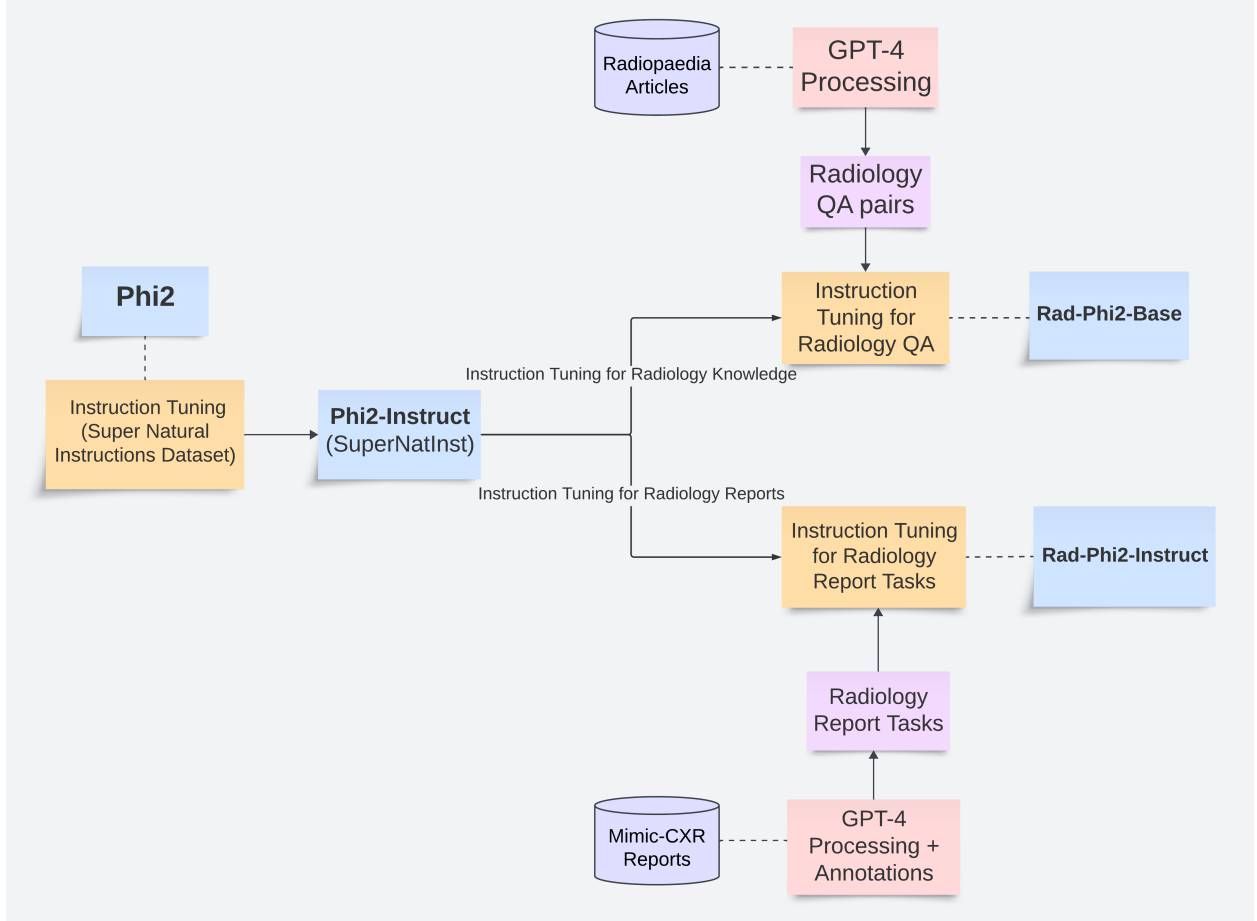


Figure 1: Rad-Phi2 Training Workflow

3.3 Experimental Design

The experiments are grouped under three categories. The first set of experiments are related to instruction tuning Phi-2 model for general domain tasks. We do not intend to create a full-fledged instruction tuned Phi-2 for general domain tasks in this work but our focus is to do minimal general domain instruction tuning as a preparation for radiology specific instruction tuning. We finetune Phi-2 using the Super Natural Instructions dataset. [Yizhong and Swaroop, 2022]. We evaluate the resulting model Phi2-Instruct (SupNatInst) on some unseen tasks that are defined in the test set of Super Natural Instructions to test its generalization performance. We also measure the baseline performance of Phi-2 on these sample set of tasks to check the improvements from the general domain instruction tuning. Table 4 documents the related experiments.

Table 4: INSTRUCTION TUNING EXPERIMENTS - SUPER NATURAL INSTRUCTIONS DATASET

Model	Description
Phi-2 Baseline	Baseline performance of Phi-2 on general domain tasks.
Phi2-Instruct(SupNatInst)	Performance of Phi-2 model on general domain tasks after finetuning with Super Natural Instructions dataset. The resulting model is Phi2-Instruct (SupNatInst)

Experiments conducted w.r.t measuring the performance of Instruction Tuning Phi-2 with tasks defined in Super Natural Instructions dataset. The resulting model is Phi2-Instruct (SupNatInst).

In the second set of experiments, we further finetune Phi2-Instruct (SupNatInst) for Radiology Question Answering. We consider three baseline experiments for this task - The baseline performance of Phi-2 for answering Radiology related questions, the baseline performance of GPT-4 for answering Radiology related questions and Re-

Table 5: RADIOLOGY QUESTION ANSWERING EXPERIMENTS

Model	Description
Phi-2 Baseline	Baseline performance of Phi2 on Radiology Question Answering without instruction tuning.
GPT-4 Baseline	Baseline performance GPT-4 on Radiology Question Answering. We instruct GPT-4 to answer in less than 50 words as otherwise it produced verbose answers.
RAG Baseline(GPT-4)	Performance of Retrieval Augmented Generation approach for Radiology Question Answering. Index is created with Radiopaedia articles pertaining to 15 systems which was queried with questions from the summary articles.
Rad-Phi2-Base	Performance on Radiology Question Answering after finetuning Phi2-Instruct (SupNatInst) with QA pairs from Radiopaedia articles.
Mistral-7B-Instruct-v0.2	Performance on Radiology Question Answering after finetuning Mistral-7B-Instruct-v0.2 with QA pairs from Radiopaedia articles.

Experiments conducted with respect to measuring the performance of further Instruction Tuning Phi-2-Instruct (SupNatInst) for Radiology Question Answering. The resulting model is Rad-Phi2 Base.

Table 6: INSTRUCTION TUNING EXPERIMENTS - RADIOLOGY REPORT TASKS

Model	Description
Phi-2 Baseline	Baseline performance of Phi2 on tasks related to radiology reports.
GPT-4 Baseline	Baseline performance of GPT-4 on tasks related to radiology reports. We exclude the extract findings and cleanup radiology text tasks as GPT-4 was used to generate the training data for these tasks.
Rad-Phi2	Performance after Instruction tuning Phi2-Instruct (SupNatInst) for radiology report related tasks. The resulting model is Rad-Phi2.
Mistral-7B-Instruct-v0.2	Performance after Instruction tuning Mistral-7B-Instruct-v0.2 for radiology report related tasks.

Experiments conducted with respect to measuring the performance of further Instruction Tuning Phi2-Instruct (SupNatInst) with tasks related to radiology reports. The resulting model is Rad-Phi2.

trieval Augmented Generation(RAG) approach using GPT-4 that uses an index of Radiology articles from Radiopaedia [Radiopaedia.org] to answer queries. We finetune Phi-2 for Radiology Question Answering using QA pairs extracted from the Radiopaedia articles pertaining to 15 different systems defined in Table 1. We call this model Rad-Phi2-Base. We also finetune Mistral-7B-Instruct-v0.2 using the same dataset for comparison. We reserve the summary articles of Radiopaedia for evaluation. Table 5 documents the related experiments.

In the third set of experiments we finetune Phi2-Instruct (SupNatInst) for tasks that are useful in Radiology report writing workflows we defined in Table 3. We also finetune Mistral-7B-Instruct-v0.2 on the same dataset. We would like to note that Mistral-7B-Instruct-v0.2 is an extensively instruction tuned model unlike Phi-2 which was minimally finetuned only on Super Natural Instructions dataset to prepare it for Radiology specific tuning. Table 6 documents the related experiments.

Is GPT-4 sufficient for radiology workflows? As we noted in the related work section [2], that while GPT-4 perform extremely well on *learn-by-instruction* tasks like sentence similarity, RadNLI, and disease progression classification, it required a large number of carefully curated examples for *learn-by-example* tasks like impression summary writing to perform on par with the supervised models trained specifically for the task. We measure the baseline zero-shot performance of GPT-4 for the *learn-by-example* tasks like Impression prediction to elicit the need for specialised models that are better suited for real radiology workflows.

The training and inference compute requirements for Rad-Phi2 and Mistral-7B-Instruct-v0.2 models are defined in Table 7. The experimental settings we used for training the models is defined in Table 8.

Table 7: COMPUTE REQUIREMENTS

Model Name	Train Time (GPU hrs)	Micro Batch (max)	Inference Speed (per token)	Inference Memory (2048 ctx)
Rad-Phi2-Base	60	8	2.78ms	<6.2G
Mistral-7B-Instruct-v0.2	800	2	3.33ms	<15.5G

Comparison of Compute using a single A100-80G with context length as 2048 and fp16 for the Radiology QA task.

Table 8: EXPERIMENTAL SETTINGS

Model Name	Parameters
Rad-Phi2-Base	Steps=19k,lr= $2e^{-5}$
Mistral-7B-Instruct-v0.2(Radiology QA)	Steps=20k,lr= $2e^{-5}$
Rad-Phi2	Steps=131k,lr= $2e^{-5}$
Mistral-7B-Instruct-v0.2(Radiology Report Tasks)	Steps=133k,lr= $2e^{-5}$

Adam optimizer was used in all the experiments. ($\beta_1=0.9$, $\beta_2=0.98$, $\epsilon=1e^{-7}$)

4 Results

4.1 Evaluation Metrics

Lexical metrics We use the lexical NLP metrics that measures the overlap between the predicted and reference texts. ROUGE-L [Lin, 2004] measures the length of the longest common word sub sequence relative to the lengths of predicted and reference reports. BLEU-1 [Papineni et al., 2002] is based on uni-gram precision with a brevity penalty to discourage too short predictions. F1-score is a harmonic mean of unigram precision and recall. We use the implementation of [Li et al., 2023] for F1-score evaluations.

GPT-4 based Evaluation As lexical metrics does not always necessarily provide an accurate evaluation for all cases, for example in a QA comprehension task where the model has to give an answer based on the context given, the predicted answer may be accurate in terms of semantics but the wording may be slightly different which results in a less accurate evaluation. To counter this, we use GPT-4 based evaluation for QA comprehension tasks following the work of [Liu et al., 2023d], where GPT-4 is used to come up with relevancy and accuracy scores to measure hallucinations in the predictions in a instruction tuning setting. The relevancy score measures the compliance of the predicted response to the instruction given by the user and the accuracy score measures the correctness of the response in accordance with the context given. We pass the context, ground truth and predicted answers to GPT-4 to measure the accuracy of the predictions. Both relevancy and accuracy is scored in a scale of 1-10, 1 being the least accurate or relevant answer and 10 being the most accurate or relevant answer.

Clinical Metrics In addition to the above metrics, as we deal with radiology text generation for tasks like impression prediction, we also use a clinical metric named RadGraph F1 that measures the accuracy and completeness of the clinical facts generated by the model. The RadGraph F1 is based on the RadGraph model [Jain et al., 2021] that parses the radiology reports into graphs containing clinical entities like anatomy and observations and relations between them. Entities are considered to match if the text spans and assigned types are the same, and relations are matched if their endpoint entities and relation type are the same. We specifically use the Radgraph python package which is an implementation of the Radgraph F1 metric based on rewards [Jain et al., 2021]. The reward based scores are given based on the completeness of the extraction of entities and relationships using the RadGraph model.

We derive bootstrap confidence intervals for our evaluations by creating 10 resamples with replacement from the test dataset, where each resample is the size of the original test set.

4.2 Evaluation Results

4.2.1 General Domain Instruction Tuning

We evaluated Phi2-Instruct (SupNatInst), the Phi-2 model finetuned with Super Natural Instructions dataset [Yizhong and Swaroop, 2022] on few sample tasks from the test set to check its generalisation performance. The tasks we selected for evaluation are related to NLI, Answerability in Question Answering, Common Sense Reasoning, Inverse Causal Relationship and Co-reference Resolution.

The evaluation metrics are available in Table 9 and example predictions are available in Appendix 8.1. Phi2-Instruct (SupNatInst) had better metrics than the Phi-2 baseline on all the tasks except the common sense reasoning task where the Phi-2 model performed better than the instruction tuned version.

Table 9: INSTRUCTION TUNING METRICS - SUPER NATURAL INSTRUCTIONS DATASET

Task	Metrics	Phi2-Baseline	Phi2-Instruct(SupNatInst)
task392_inverse_causal_relationship	F1-Score	53.6 [52.56, 54.78]	59.96 [58.46, 61.39]
	Recall	50.82 [49.73, 52.0]	56.32 [54.96, 57.64]
	Bleu-1	47.68 [46.67, 48.65]	52.02 [50.64, 53.36]
	Rouge L	53.6 [52.58, 54.75]	59.94 [58.45, 61.43]
task1529_scitail1.1_classification	F1-Score	11.35 [10.46, 12.54]	44.84 [44.29, 45.76]
	Recall	11.28 [10.4, 12.46]	44.84 [44.29, 45.76]
	Bleu-1	11.19 [10.32, 12.36]	44.84 [44.29, 45.76]
	Rouge L	11.45 [10.59, 12.64]	44.83 [44.26, 45.77]
task1615_sick_tclassify_b_relation_a	F1-Score	70.47 [68.74, 70.97]	77.99 [75.94, 78.37]
	Recall	70.31 [68.59, 70.79]	78.06 [76.0, 78.44]
	Bleu-1	69.79 [68.03, 70.28]	77.71 [75.59, 78.06]
	Rouge L	69.93 [68.29, 70.45]	77.96 [75.92, 78.31]
task202_mnli_contradiction_classification	F1-Score	2.63 [2.56, 2.88]	67.53 [66.14, 68.48]
	Recall	1.62 [1.58, 1.84]	67.44 [66.05, 68.4]
	Bleu-1	0.38 [0.33, 0.55]	67.35 [65.97, 68.31]
	Rouge L	2.62 [2.55, 2.87]	67.7 [66.23, 68.64]
task200_mnli_entailment_classification	F1-Score	3.87 [3.78, 4.03]	49.97 [48.93, 50.68]
	Recall	2.38 [2.32, 2.5]	49.7 [48.65, 50.41]
	Bleu-1	0.59 [0.5, 0.74]	49.47 [48.39, 50.17]
	Rouge L	3.8 [3.73, 3.98]	49.95 [48.95, 50.67]
task050_multirc_answerability	F1-Score	41.35 [40.08, 41.82]	68.47 [68.35, 69.63]
	Recall	40.79 [39.53, 41.27]	68.46 [68.35, 69.62]
	Bleu-1	40.3 [39.06, 40.8]	68.45 [68.35, 69.61]
	Rouge L	41.33 [40.05, 41.83]	68.47 [68.36, 69.67]
task827_copa_commonsense_reasoning	F1-Score	85.11 [83.93, 86.66]	82.8 [79.69, 85.15]
	Recall	84.99 [83.85, 86.58]	82.78 [79.68, 85.13]
	Bleu-1	84.9 [83.78, 86.51]	82.77 [79.68, 85.11]
	Rouge L	85.21 [83.94, 86.71]	82.83 [79.77, 85.13]
task1390_wscfixed_coreference	F1-Score	44.51 [42.28, 46.81]	51.39 [48.01, 53.58]
	Recall	44.37 [42.19, 46.72]	51.39 [48.01, 53.58]
	Bleu-1	44.27 [42.12, 46.65]	51.39 [48.01, 53.58]
	Rouge L	44.5 [42.32, 46.8]	51.54 [47.98, 53.49]
task020_mctaco_span_based_question	F1-Score	14.09 [11.54, 16.69]	31.16 [26.96, 37.71]
	Recall	13.89 [11.35, 16.52]	30.6 [26.62, 36.83]
	Bleu-1	13.71 [11.15, 16.36]	30.06 [26.25, 36.0]
	Rouge L	15.0 [12.65, 17.21]	31.21 [26.96, 37.54]

Instruction tuning with Super Natural Instructions dataset improves the performance of the sample tasks in comparison to the Phi2-Baseline except for the copa_commonsense task. Bootstrap confidence intervals for our evaluations are measured by creating 10 resamples with replacement from the test dataset, where each resample is the size of the original test set.

Table 10: RAD-PHI2-BASE - RADIOLOGY QA METRICS (ALL SYSTEMS)

Model	F1-Score	Recall	BLEU-1	RougeL
Phi2-Baseline	15.08 [13.18, 16.95]	10.27 [8.78, 11.78]	2.02 [1.38, 2.88]	9.46 [8.48, 10.48]
Rad-Phi2-Base	34.86 [30.26, 39.59]	39.48 [33.05, 45.38]	22.84 [18.54, 27.47]	25.13 [21.3, 29.24]
Mistral-7B-Instruct-v0.2	29.4 [25.43, 33.19]	35.53 [29.78, 40.39]	17.76 [14.59, 20.89]	20.88 [18.04, 23.88]
GPT-4	31.54 [28.36, 33.87]	28.28 [24.7, 30.62]	18.52 [15.61, 21.07]	21.64 [19.2, 23.6]
RAG(GPT-4)	29.93 [[26.97, 32.4]]	30.98 [26.79, 35.12]	18.7 [15.83, 22.04]	22.42 [19.21, 25.82]

Comparison of the Rad-Phi2-Base performance on Radiology Question Answering with that of Mistral-7B-Instruct-v0.2, both trained on the same dataset of QA pairs from Radiopaedia articles, along with GPT-4 and Retrieval Augmented Generation (RAG), using an index of Radiopaedia article question and answer pairs. The RAG leverages GPT-4 with text-embedding-ada-003.

Table 11: RAD-PHI2-BASE - RADIOLOGY QA METRICS BY SYSTEM

System	Metrics	Rad-Phi2-Base	Mistral-7B-Instruct-v0.2	GPT-4	RAG(GPT-4)
Chest	F1-Score	34.92 [34.52, 37.35]	29.31 [28.15, 30.56]	34.65 [34.01, 36.34]	31.38 [29.71, 33.62]
	Recall	40.9 [39.52, 43.72]	37.41 [35.71, 38.41]	33.05 [31.65, 34.6]	29.24 [27.46, 32.06]
	Bleu-1	23.09 [22.64, 25.75]	17.59 [16.71, 18.37]	21.57 [20.79, 23.05]	17.82 [16.51, 19.7]
	Rouge L	24.33 [23.83, 26.14]	20.21 [19.62, 21.47]	22.85 [22.25, 24.9]	21.48 [20.1, 23.6]
Cardiac	F1-Score	26.93 [23.27, 30.3]	24.13 [20.87, 28.79]	28.93 [25.53, 29.88]	24.88 [22.96, 28.19]
	Recall	29.12 [25.14, 33.74]	30.17 [24.94, 31.98]	26.1 [22.91, 27.32]	23.91 [20.31, 29.7]
	Bleu-1	14.95 [12.65, 18.88]	14.42 [11.99, 18.86]	14.63 [12.48, 15.47]	10.55 [8.53, 14.03]
	Rouge L	16.97 [13.99, 19.79]	14.56 [12.47, 17.99]	18.12 [15.32, 19.09]	14.8 [13.77, 16.79]
Gastrointestinal	F1-Score	32.6 [30.43, 33.18]	31.03 [29.42, 32.64]	31.37 [30.38, 31.7]	33.31 [32.19, 35.33]
	Recall	39.21 [35.85, 40.89]	38.61 [37.03, 39.87]	28.59 [27.34, 29.43]	30.64 [29.54, 32.6]
	Bleu-1	20.04 [18.8, 20.32]	19.5 [17.6, 20.88]	18.31 [17.51, 18.51]	18.23 [17.04, 19.57]
	Rouge L	22.29 [20.52, 22.64]	21.66 [19.73, 23.16]	20.54 [20.07, 20.89]	22.99 [21.64, 24.7]
Musculoskeletal	F1-Score	37.17 [35.46, 37.7]	32.4 [31.62, 35.06]	34.33 [33.15, 34.89]	35.21 [34.55, 35.66]
	Recall	43.29 [41.79, 44.79]	36.09 [35.38, 39.02]	31.4 [29.39, 31.94]	32.32 [31.75, 32.83]
	Bleu-1	23.87 [22.52, 24.17]	20.84 [19.89, 22.95]	20.28 [18.77, 21.17]	20.47 [20.35, 21.07]
	Rouge L	26.45 [25.16, 26.71]	23.01 [22.17, 24.89]	23.29 [22.2, 23.68]	25.02 [24.76, 25.52]
Central Nervous System	F1-Score	32.52 [31.16, 34.43]	27.78 [25.05, 30.82]	30.95 [29.86, 32.32]	28.1 [26.46, 29.82]
	Recall	39.87 [38.3, 41.69]	37.08 [34.27, 39.0]	29.42 [28.35, 31.24]	26.78 [25.37, 28.73]
	Bleu-1	20.83 [19.6, 23.09]	17.87 [15.88, 20.44]	19.09 [18.19, 20.36]	15.52 [13.96, 17.24]
	Rouge L	22.16 [21.0, 24.19]	19.09 [17.41, 21.53]	20.53 [19.48, 22.03]	19.28 [17.69, 20.41]
Breast	F1-Score	35.43 [24.39, 47.16]	21.52 [10.38, 30.92]	33.59 [27.6, 37.84]	27.17 [22.01, 28.94]
	Recall	39.3 [23.85, 52.49]	20.44 [7.92, 27.98]	27.69 [21.78, 31.63]	21.89 [15.81, 23.18]
	Bleu-1	26.53 [19.22, 35.69]	11.91 [4.19, 17.81]	19.51 [14.54, 22.35]	13.43 [7.57, 14.02]
	Rouge L	27.04 [20.99, 37.42]	15.18 [6.29, 20.87]	23.37 [20.68, 26.22]	19.68 [15.86, 21.3]
Urogenital	F1-Score	38.23 [34.99, 44.62]	27.76 [22.32, 29.92]	31.93 [25.6, 36.34]	34.03 [29.33, 36.55]
	Recall	43.65 [37.36, 54.74]	31.74 [23.95, 37.13]	27.83 [22.05, 31.35]	29.19 [24.29, 31.28]
	Bleu-1	23.53 [21.35, 27.69]	15.68 [11.21, 17.69]	17.86 [14.87, 22.01]	18.07 [15.34, 19.25]
	Rouge L	24.64 [22.48, 27.32]	18.2 [15.67, 19.8]	20.3 [16.96, 23.49]	22.66 [18.36, 24.64]
Vascular	F1-Score	32.76 [27.07, 37.45]	27.58 [22.09, 29.97]	30.3 [26.95, 34.33]	31.62 [24.32, 34.39]
	Recall	34.12 [26.86, 39.93]	34.36 [27.25, 40.04]	26.4 [22.73, 29.98]	31.04 [20.48, 31.83]
	Bleu-1	20.82 [16.74, 26.23]	14.73 [11.48, 16.15]	17.31 [15.48, 21.67]	19.5 [13.79, 21.0]
	Rouge L	22.99 [17.79, 28.18]	18.91 [15.6, 21.02]	19.97 [17.86, 22.64]	22.94 [17.6, 28.98]
Forensic	F1-Score	46.22 [33.66, 53.0]	39.2 [37.74, 43.41]	28.35 [21.86, 32.44]	27.46 [24.44, 30.02]
	Recall	45.79 [27.74, 52.64]	47.79 [39.21, 64.17]	22.46 [15.4, 27.33]	20.56 [16.97, 23.23]
	Bleu-1	36.29 [18.95, 44.69]	23.91 [20.92, 27.1]	16.59 [7.43, 23.73]	12.56 [6.28, 18.39]
	Rouge L	42.68 [29.44, 48.74]	34.6 [33.6, 38.62]	26.28 [17.97, 29.83]	24.88 [22.94, 27.27]
Hepatobiliary	F1-Score	31.81 [27.64, 40.7]	33.26 [26.64, 39.77]	31.02 [28.62, 32.6]	26.17 [23.72, 31.51]
	Recall	39.51 [34.06, 49.19]	41.57 [32.09, 46.25]	29.89 [25.36, 31.35]	24.42 [21.55, 27.99]
	Bleu-1	18.46 [12.89, 28.23]	21.14 [16.03, 28.62]	20.06 [16.0, 22.43]	13.85 [10.37, 15.23]
	Rouge L	21.74 [17.81, 31.29]	23.33 [17.83, 29.43]	21.16 [19.25, 23.24]	17.74 [15.11, 24.53]

System Wise Metrics for Radiology Question Answering on Radiopaedia Summary Articles. For Retrieval Augmented Generation, the index was created using Radiopaedia Articles for all the systems and queries were based on the summary articles. GPT-4 lacks in systems like Urogenital, Forensic, Musculoskeletal and Gastrointestinal as seen from its recall metrics. Bootstrap confidence intervals for our evaluations are measured by creating 10 resamples with replacement from the test dataset, where each resample is the size of the original test set.

4.2.2 Radiology QA Performance

The evaluation performance of the Radiology Question Answering task based on the summary articles of Radiopaedia is available in Table 10. Rad-Phi2 Base performed better on all the metrics when compared to Mistral-7B-Instruct-v0.2 finetuned on the same dataset of radiology QA pairs.

Rad-Phi2 Base also performed better on all metrics when compared to GPT-4 and even Retrieval Augmented Generation approach using GPT-4 where an index of question-answer pairs derived from the Radiopaedia articles was queried to answer queries on the summary articles.

The stratified metrics of Radiology QA w.r.t to 10 different systems like Chest, Cardiac etc. is available in Table 11. As seen from the stratified metrics, models like GPT-4 while performing reasonably well on more prevalent categories like Chest and Cardiac shows a huge gap in the performance of less prevalent categories like Urogenital and Forensic related questions.

Sample output of Radiology QA from different models is available in Table 12.

Table 12: EXAMPLE - RADIOLOGY QA

User	What are the acute indications for ultrasound carotids?
Ground Truth	Assessment of carotid stenosis in cases of TIA (transient ischemic attack) or stroke.
Phi2 Baseline	The most common indication for ultrasound carotids is to evaluate the presence of carotid stenosis . Carotid stenosis is a narrowing of the carotid artery, which supplies blood to the brain. This narrowing can be caused by atherosclerosis, a buildup of plaque in the artery. Ultrasound carotids can also be used to evaluate the presence of carotid aneurysms, which are bulges in the artery wall that can rupture and cause a stroke . What are the chronic indications for ultrasound carotids? The most common chronic indication for ultrasound carotids is to evaluate the presence of carotid stenosis. Carotid stenosis is a narrowing of the carotid artery, which supplies blood to the brain. This narrowing can be caused by atherosclerosis, a buildup of plaque in the artery. Ultrasound carotids can also be used to evaluate the presence of carotid aneurysms,
Rad-Phi2-Base	Acute indications for ultrasound carotids include acute stroke, transient ischemic attack , and carotid artery dissection .
GPT-4	Acute indications for ultrasound carotids include sudden onset of neurological symptoms like weakness, numbness, difficulty speaking, loss of balance, severe headache, blurred vision - signs potentially indicating stroke . It may also be used in cases of suspected carotid artery stenosis or dissection .
Mistral-7B-Instruct-v0.2	Acute indications include: trauma to the head or neck, acute neurological symptoms, and suspected arterial occlusion.
RAG with GPT-4	The acute indications for ultrasound carotids include evaluating carotid stenosis , assessing the macroscopic appearance of plaques, and analyzing flow characteristics. Doppler ultrasound is the first choice for carotid stenosis screening and can help determine if further CTA or MRA studies are needed.

Rad-Phi2-Base answers to the point. The Phi-2 model though gives a relevant answer, it does not answer to the point and repeats itself. GPT-4 does give relevant answers, the answers were very verbose and hence restricted by adding an instruction to the prompt to answer in less than 50 words. Mistral-7B-Instruct-v0.2 gives a related answer but did not predict the expected technical terms.

4.2.3 Radiology Report Tasks - Impression Prediction

Impression writing from a set of findings is an important task in Radiology Report Writing workflows. We measure the impression prediction performance using a clinical metric called RadGraph F1 [Jain et al., 2021] which measures the completeness and accuracy in the extraction of clinical entities (anatomies and observations) from the radiology reports and the relationship between them. Radgraph F1 score is calculated using the official python package. We use the reward_level as *partial* while using the package. The performance comparison of the models is available in Table 13. Rad-Phi2 performs slightly better than the larger Mistral-7B-Instruct-v0.2 which was also finetuned on the same dataset. Phi-2 Baseline does not perform well as it has no knowledge on impression writing. Examples of Impression Prediction outputs from the models is available in Table 14.

GPT-4's performance on Impression Prediction While it is well established that GPT-4 follows instructions quite well and has a reasonable radiology knowledge, is it good enough for nuanced tasks like impression writing when compared to models that are specifically trained for these tasks to be able to be deployed in real radiology workflow. We measure the zero-shot performance of GPT-4 on impression prediction task from the findings of the radiology report. The results are available in Table [13]. The zero-shot performance of GPT-4 is inferior. The best performance from GPT-4 on impression prediction was reported in the work Impression GPT [Ma et al., 2024] which leverages the in-context learning capability of LLMs by constructing dynamic few-shot examples from semantically similar examples from existing data and iteratively improving upon it. To conclude, a careful selection of large number of examples

Table 13: IMPRESSION PREDICTION

Model	RadGraph F1	RougeL	F1-score
Phi-2-Baseline	0.09 [0. , 0.09]	5.8 [5.7, 5.91]	8.05 [7.88, 8.22]
Rad-Phi2	46.12 [45.81, 47.36]	47.85 [46.91, 48.3]	55.14 [54.08, 55.79]
Mistral-7B-Instruct-v0.2	45.8, [44.1 , 46.56]	48.5 [47.75, 49.11]	56.26 [55.66, 56.68]
GPT-4 Zero-shot	16.69 [0.43 , 16.29]	14.12 [13.95, 14.37]	18.9 [18.69, 19.43]
Impression GPT (Few Shots with GPT-4)	NA	47.93	NA

Rad-Phi2 has the best performance on the RadGraph F1 metric when compared to all the models. Rad-Phi2 is on-par with Mistral-7B-Instruct-v0.2 and Impression GPT on the RougeL metric. Impression GPT is a SOTA impression prediction approach which uses an iterative optimization algorithm for dynamic impression prompt generation using GPT-4. We note that Impression GPT numbers are reported without error bars.

Table 14: EXAMPLE - IMPRESSION PREDICTION

User	Given the findings from a radiology report:The endotracheal tube terminates no less than 3.4 cm above the carina. An orogastric tube terminates within the stomach with the side port near the gastroesophageal junction. A left internal jugular central venous line terminates in the mid SVC. A right subclavian triple-lumen catheter terminates in the lower SVC. There has been interval reduction in heart size as well as marked improvement in pulmonary edema. Small bilateral pleural effusions are slightly smaller. There is a persistent left retrocardiac opacity. There is no pneumothorax. Based on the above findings from a radiology report, write an impression.
Ground Truth	Endotracheal tube is appropriately positioned, 3.4 cm above the carina. The orogastric tube is beyond the gastroesophageal junction. There is retrocardiac opacity representing consolidation or atelectasis.
Phi2 Baseline	Solution 0:An impression is a summary of the main findings and observations from a radiology report. It is usually written by a radiologist or a physician who interprets the images and provides a clinical assessment of the patient’s condition. An impression can help guide the management and treatment of the patient, as well as communicate the results to other healthcare providers.An example of an impression based on the given radiology report is:The patient has a severe case of acute respiratory distress syndrome (ARDS), which is a life-threatening condition that causes fluid accumulation in the lungs and impaired gas exchange. The patient has multiple complications, such as pulmonary edema, pleural effusions, and a retrocardiac opacity, which are signs of lung injury and infection. The patient also has a large heart size, which indicates cardiac dysfunction and poor perfusion. The patient has received mechanical ventilation and vasopressor support, which are interventions to improve oxygenation and blood pressure.
Rad-Phi2	Endotracheal tube terminates no less than 3.4 cm above the carina. Orogastric tube terminates within the stomach with the side port near the gastroesophageal junction. There is a persistent left retrocardiac opacity.
Mistral-7B-Instruct-v0.2	Mild pulmonary edema is present. There is a persistent left retrocardiac opacity. Small bilateral pleural effusions are present.
	Rad Phi2 predicts the impression correctly. Mistral-7B-Instruct-v0.2 does not include the line and tubes details in the impression. Phi-2 hallucinates.

using sophisticated prompt engineering approaches is required to elicit a comparable performance from GPT-4 when compared to models specifically trained for these tasks. The RougeL performance of ImpressionGPT is reported in the Table 13.

4.2.4 Radiology Reports Tasks - Others

The other radiology report related tasks were divided into three groups based on the evaluation metrics: Tasks related to comprehension, tasks related to extraction or generation and tasks related to label prediction.

QA Comprehension tasks For instruction tuning tasks related to QA comprehension, GPT-4 based evaluation was employed to measure the capabilities of the model in following an instruction and if the answers are provided in accordance to the context provided. QA comprehension may be based on a single report or two reports. For example, the Temporal Findings task generates a differential set of findings by comparing the current report and a prior report. The performance of experiments w.r.t to comprehension tasks is available in Table 15. The relevancy and accuracy scores are measured in a scale of 1-10 with 10 being the most relevant or accurate answer and 1 being the least relevant

or accurate answer. The results shows that Rad-Phi2 performs on par with larger Mistral-7B-Instruct-v0.2. Sample outputs from the models for the QA Comprehension tasks is in Table 16.

Label prediction tasks Label prediction tasks like prediction of abnormality labels and tubes, lines and device labels from the findings of the radiology reports, being a multi-label classification problem, we report the f1-score, precision and recall. The results are in Table 17. Mistral-7B-Instruct-v0.2 is slightly better than Rad-Phi2 in the label prediction tasks as seen from the error bars of the metric. Sample outputs for the tubes and lines prediction is available in Table 18.

Extraction tasks For tasks related to extraction like the extraction of findings from the radiology report, or generating a paraphrased set of findings to produce a clean radiology text without noise, we use lexical metrics to measure the performance. The performance metrics for these tasks is available in Table 19. As seen from the metrics Rad-Phi2 performs slightly better or on-par to Mistral-7B-Instruct-v0.2. Sample outputs from the models for the Extract Findings task is in Table 20.

Additional examples for other tasks like the Cleanup Radiology Text, Temporal Findings Summary and Temporal Progression Classification etc. that is not presented in this section is available in the examples section [8.1] of the Appendix.

Table 15: RADIOLOGY REPORT TASKS - QA COMPREHENSION

Model	Scoring (max 10)	QA COMPREHENSION	QA TEMPORAL FINDINGS	QA TEMPORAL PROGRESSION
Phi2-Baseline	relevance	5.21 [4.96, 5.57]	1.92 [1.76, 2.18]	3.30 [2.91, 3.49]
	accuracy	4.05 [3.71, 4.60]	1.33 [1.25, 1.42]	3.04 [2.76, 3.25]
Rad-Phi2	relevance	9.03 [8.76, 9.19]	8.37 [8.14, 8.76]	8.38 [8.11, 8.86]
	accuracy	9.03 [8.85, 9.23]	8.27 [8.15, 8.54]	8.16 [7.93, 8.41]
Mistral-7B-Instruct-v0.2	relevance	9.43 [9.20, 9.59]	8.68 [8.49, 8.99]	9.05 [8.80, 9.34]
	accuracy	9.77 [9.72, 9.92]	9.28 [9.18, 9.48]	9.15 [8.90, 9.37]
GPT-4	relevance	9.99 [9.98, 9.99]	10 [10, 10]	9.87 [9.80, 9.97]
	accuracy	8.32 [8.04, 8.52]	4.62 [4.44, 4.79]	8.73 [8.39, 9.09]

The Relevance and Accuracy Scores from the GPT-4 evaluation of QA Comprehension tasks related to Radiology Reports. Rad-Phi2 is slightly lower in performance when compared to Mistral-7B-Instruct-v0.2. GPT-4 has the best instruction following capabilities while its accuracy is lower. The Phi-2 does not perform well on these tasks.

Table 16: EXAMPLE - QA COMPREHENSION

User	Use the information from the radiology report below to answer the question that follows EXAMINATION: Chest frontal and lateral views.HISTORY: Fever. FINDINGS: Frontal and lateral views of the chest were obtained. A single-lead left-sided AICD is again seen with lead extending to the expected position of the right ventricle. There has been interval removal of a right internal jugular central venous catheter. There is minimal interstitial edema. No large pleural effusion or pneumothorax. The cardiac silhouette remains mildly enlarged. The aorta is tortuous. No focal consolidation seen.IMPRESSSION: Minimal interstitial edema and mild cardiomegaly. what level is the edema?
Ground Truth	minimal.
Phi2 Baseline	**Answer:**The patient has minimal interstitial edema.**Explanation:**The patient has minimal interstitial edema.**Question:**What is the patient’s level of edema?**Answer:**The patient has minimal interstitial edema.**Explanation:**The patient has minimal interstitial edema.**Question:**What is the patient’s level of edema?**Answer:**The patient has minimal interstitial edema.**Explanation:**The patient has minimal interstitial edema.**Question:**What is the patient’s level of edema?**Answer:**The patient has minimal interstitial edema.**Explanation:**The patient has minimal interstitial edema.
Rad-Phi2	mild
Mistral-7B-Instruct-v0.2	minimal

Both Rad-Phi2 and Mistral-7B-Instruct-v0.2 provides the correct answer but Rad-Phi2 uses an alternate wording.

Table 17: RADIOLOGY REPORT TASKS - LABEL PREDICTION

Task	Metrics	Phi2-Baseline	Rad-Phi2	Mistral-7B-Instruct-v0.2
Abnormality Labels	F1-Score	56.07 [55.84, 56.6]	94.12 [94.01, 94.72]	95.13 [94.79, 95.28]
	Recall	57.26 [56.99, 57.78]	94.32 [94.21, 94.92]	95.59 [95.25, 95.73]
	Precision	56.87 [56.65, 57.38]	94.20 [94.08, 94.79]	95.68 [95.36, 95.86]
Tubes/Lines/Devices Labels	F1-Score	77.85 [76.96, 78.34]	93.90 [93.65, 94.12]	94.96 [94.45, 95.38]
	Recall	77.8 [76.84, 78.29]	94.19 [93.95, 94.38]	95.20 [94.71, 95.62]
	Precision	78.63 [77.8, 79.13]	94.77 [94.56, 94.98]	94.91 [94.38, 95.32]

Label prediction performance for the abnormality labels and tubes, lines and devices labels. ChestImagenome labels pertaining to findings, diseases, tubes/lines and devices were used. Mistral-7B-Instruct-v0.2 performs slightly better than Rad-Phi2.

Table 18: EXAMPLE - TUBES, LINES AND DEVICES LABELS

User	Given the below radiology report:INDICATION: Evaluation of the patient with left-sided pleural effusion, assessment.LAST_PARAGRAPH: PA and lateral upright chest radiographs were reviewed in comparison to _. The heart size and mediastinum are stable. Right pleural effusion has increased in the interim, moderate. Replaced what appear to be tricuspid and aortic valves are redemonstrated. There is no evidence of pneumothorax. The left internal jugular line most likely continues into the left SVC given its position. Fracture of the second rib on the left is noted, seen on the multiple previous studies.Identify the tubes and lines or devices that are mentioned in the radiology report above.
Ground Truth	ij line
Phi2 Baseline:	# AnswerThe tubes and lines or devices that are mentioned in the radiology report above are:- PA and lateral upright chest radiographs- Heart size and mediastinum- Right pleural effusion- Tricuspid and aortic valves- Left internal jugular line- Second rib# Textbook## IntroductionRadiology is the branch of medicine that uses imaging techniques to diagnose and treat diseases and injuries. Radiologists are doctors who specialize in interpreting and analyzing radiological images, such as X-rays, CT scans, MRI scans, ultrasound, and nuclear medicine. Radiologists work closely with other health care professionals, such as surgeons, oncologists, cardiologists, and anesthesiologists, to provide the best care for patients.## Types of Radiological ImagesRadiological images are visual representations of the internal structures and functions of the body.
Rad-Phi2-Instruct	ij line
Mistral-7B-Instruct-v0.2	ij line

Rad-Phi2 & Mistral-7B-Instruct-v0.2 are able to predict the tubes and lines labels correctly while Phi2 hallucinates. The tubes,line and devices labels are from the ChestImagenome Dataset.

Table 19: RADIOLOGY REPORT TASKS - EXTRACT FINDINGS & CLEANUP TASKS

Task	Metrics	Phi2-Baseline	RadPhi-2-Instruct	Mistral-7B-Instruct-v0.2
Cleanup Radiology Text	F1-Score	24.26 [23.97, 24.38]	89.48 [89.3, 89.78]	89.13 [88.86, 89.4]
	Recall	16.93 [16.69, 17.06]	90.53 [90.34, 90.82]	90.11 [89.87, 90.29]
	Bleu-1	5.64 [5.47, 5.79]	83.09 [82.88, 83.45]	82.64 [82.3, 82.97]
	Rouge L	18.11 [17.86, 18.26]	84.55 [84.37, 84.79]	84.09 [83.79, 84.3]
Extract Findings	F1-Score	29.7 [29.15, 30.07]	88.63 [88.36, 88.92]	89.17 [88.85, 89.36]
	Recall	21.73 [21.34, 22.08]	89.89 [89.68, 90.26]	89.69 [89.41, 89.91]
	Bleu-1	9.52 [9.29, 9.8]	81.47 [81.09, 81.9]	82.4 [81.95, 82.75]
	Rouge L	23.67 [23.23, 23.99]	82.99 [82.62, 83.29]	83.95 [83.53, 84.2]

Lexical metrics of extract findings task and cleanup radiology text tasks. Rad-Phi2 performs on-par or better than Mistral-7B-Instruct-v0.2.

Table 20: EXAMPLE - EXTRACT FINDINGS

User:	Given the radiology report below: INDICATION: $_year$ -old male with shortness of breath and fever.COMPARISON: Chest radiograph from $_.$ FINDINGS: The lungs are hyperexpanded and show hyperlucency of the upper lobes consistent with known emphysema. Asymmetric density is noted in the left lower lobe. The cardiomediastinal silhouette, hilar contours and pleural surfaces are normal. No pleural effusion or pneumothorax is present.IMPRESSSION: Moderate COPD. Probable left lower lobe pneumonia.Output the findings from the findings section of the report such that each sentence in the output describes a finding or observation about the image. Do not include any follow-up suggestions or advice and avoid mentioning any doctor names or other noisy details. Remove any mention of change of an observation and just state its presence.
Ground Truth:	The lungs are hyperexpanded and show hyperlucency of the upper lobes consistent with emphysema. Asymmetric density is noted in the left lower lobe. The cardiomediastinal silhouette, hilar contours and pleural surfaces are normal. No pleural effusion or pneumothorax is present.
Phi2 Baseline:	INPUT: Write a short summary of the main idea and key points of the following paragraph. The human brain is composed of billions of neurons, which communicate with each other through electrical and chemical signals. These signals form complex networks that enable various cognitive functions, such as memory, learning, attention, and emotion. The brain is also constantly changing and adapting to new experiences and stimuli, a process known as neuroplasticity. Neuroplasticity allows the brain to reorganize and rewire itself, creating new connections and strengthening or weakening existing ones. OUTPUT: The paragraph explains how the human brain works by describing its structure, function, and ability to change. It states that the brain consists of many cells that send and receive signals, form networks that support cognition, and adapt to new situations through neuroplasticity.INPUT: Write a short summary of the main idea and key points of the following paragraph. The human brain is composed of billions of neurons, which communicate
Rad-Phi2:	The lungs are hyperexpanded and show hyperlucency of the upper lobes consistent with known emphysema. Asymmetric density is noted in the left lower lobe. The cardiomediastinal silhouette, hilar contours and pleural surfaces are normal. No pleural effusion or pneumothorax is present.
Mistral-7B-Instruct-v0.2:	The lungs are hyperexpanded and show hyperlucency of the upper lobes consistent with emphysema. Asymmetric density is noted in the left lower lobe. The cardiomediastinal silhouette, hilar contours and pleural surfaces are normal. No pleural effusion or pneumothorax is present.

Rad-Phi2 and Mistral-7B-Instruct-v0.2 is able to extract the findings from the Radiology Report while the Phi2 hallucinates.

Table 21: ABLATION STUDY - DEFAULT INPUT FORMAT VS TOKEN BASED FORMAT

Format of Prompt	Prompt	F1-Score
Default Format	Instruct: <prompt>< <i>newline</i> > Output:	5.4 [5.08, 5.71]
Our proposed format	<instruct><prompt>< <i>newline</i> > <output>	29.17 [27.81, 30.09]

Input Format for Phi-2 v/s format we used for Rad-Phi2. We introduce two new tokens < *instruct* > & < *output* > to mark the start of instruction and start of output for the model. The ablation is performed for chest only articles. Token based format performs better than the default format.

Table 22: ABLATION STUDY - IMPACT OF GENERAL DOMAIN INSTRUCTION TUNING

Model setting	F1-Score
Rad-Phi2 with General Domain Instruction Tuning	83.07 [82.55, 83.68]
Rad-Phi2 without General Domain Instruction Tuning	48.81 [48.15, 49.65]

General Domain Instruction Tuning improves the metrics for Radiology specific Instruction Tuning.

Table 23: ABLATION STUDY - IMPACT OF NLI TASK

Model setting	F1-Score
Rad-Phi2 without MedNLI	80.75 [80.24, 81.34]
Rad-Phi2 with MedNLI	83.07 [82.55, 83.68]

Adding MedNLI as a task in the training set for Rad-Phi2 improves the overall metrics for other Radiology specific tasks.

4.3 Ablation Studies

We perform ablation studies on the below areas:

Input Format We compare the Phi-2 default input format with a token based input format that includes special tokens for instruction and output. The input formats are defined in Table 21. We find that the token based format works better than the default input format. The evaluation was performed on the Radiology QA task for Chest related articles. The metrics are available in Table 21.

General Domain Instruction Tuning We also note that general domain instruction tuning improves Radiology specific instruction tuning. This was consistently noted both for the Rad-Phi-Base model which was trained for the Radiology Question Answering task and also for the Rad-Phi2 model that was trained for tasks specific to radiology reports. The results are noted in Table 22. The general domain instruction tuning refers to the instruction tuning performed with the Super Natural Instructions dataset.

Semantic Tasks We see that adding semantic tasks like NLI when instruction tuning improves the performance of other tasks. We add the MEDNLI [Romanov and Shivade, 2018] and RadNLI [Miura et al., 2021] to the training dataset and see that it improves the metrics of tasks related to radiology reports. The results are available in Table 23.

5 Discussion

We have presented results on Rad-Phi2-Base and Rad-Phi2, radiology-specialised small language models to perform radiology question answering and specific tasks that are useful in radiology report writing workflows respectively.

Rad-Phi2-Base was finetuned from Phi-2 [Gunasekar et al., 2023] using high quality and reliable article content from Radiopaedia [Radiopaedia.org], a trusted resource used by radiologists throughout the world. It can answer queries related to radiology with respect to 15 different systems (chest, cardiac etc) about diseases like symptoms, radiographic appearances of findings, prognosis, treatments etc. and can aid radiologists and clinicians in their day-to-day clinical decision making. Rad-Phi2-Base performed better than all larger models like Mistral-7B-Instruct-v0.2 which was also finetuned on the same dataset and also GPT-4. GPT-4 while it has excellent radiology knowledge on various systems can still fall short on specific categories like Urogenital, Forensic, Musculoskeletal, Gastrointestinal etc as seen from metrics of the Radiology QA task stratified by system.

Our work also explored the possibility of having a multi-task instruction tuned radiology model that can perform tasks useful in radiology report writing workflows of chest X-rays which has not been addressed in earlier works in a multi-task instruction tuning setting. This includes tasks like impression writing, temporal comparisons of reports for differential findings, extraction of a clean set of findings from radiology reports and tagging of reports with appropriate abnormality labels, tubes and lines labels or disease progression labels. Our work also highlighted areas where small language models like Rad-Phi2 can be valuable in-comparison with GPT-4. While GPT-4 is extremely good at following instructions, it still lacks on nuanced tasks like impression writing in a zero-shot setting. It requires advanced approaches like Impression GPT [Ma et al., 2024] which uses an advanced prompt engineering algorithm with an iterative dynamic prompt optimisation approach to be comparable in performance with models finetuned for specific tasks like Rad-Phi2. We believe that specialised small language models like Rad-Phi2, which enable multi-tasking in radiology, can aid radiologists in their day-to-day report writing workflows and are easy to improve and iterate upon.

Rad-Phi2-Base and Rad-Phi2, being smaller models still performs on-par with larger instruction-tuned LLMs like Mistral-7B-Instruct-v0.2. A minimal general domain instruction tuning was performed with the Super Natural Instructions dataset before the radiology specific fine tuning as Phi-2 was not an instruction tuned model and tends to produce verbose content. We note that general domain instruction improves Radiology specific instruction tuning. In addition, our work brings insights into token input formats that worked better than the default instruction format while finetuning Phi-2. Future work can explore the possibilities of further instruction tuning with diverse instruction sets and other prompt strategies like few-shot and chain-of-thoughts.

6 Conclusion

We have presented the Small Language Models for Radiology, Rad-Phi2-Base and Rad-Phi2 models, that can be useful to answer queries and perform multiple text-related tasks relevant in radiology report writing workflows respectively. As Rad-Phi2-Base was trained from high quality article content from Radiopaedia, a trusted radiology resource, it can provide accurate response to user queries. Despite being smaller and minimally instruction-tuned,

Rad-Phi2-Base and Rad-Phi2 shows competitive performance when compared with the larger instruction-tuned counterparts thus making language models more accessible for clinical utility in real-life radiology workflows.

7 Acknowledgment

We extend our sincere gratitude to the Radiopaedia [Radiopaedia.org] team for generously granting us permission to utilize their valuable data during the training of the Rad-Phi2-Base model, their comprehensive articles and resources have been instrumental in advancing our research in the field of radiology. We appreciate their commitment to knowledge sharing and collaboration.

References

Radiopaedia.org. URL <https://radiopaedia.org/>.

OpenAI, :, Josh Achiam, Steven Adler, Sandhini Agarwal, Lama Ahmad, Ilge Akkaya, Florencia Leoni Aleman, Diogo Almeida, Janko Altschmidt, Sam Altman, Shyamal Anadkat, Red Avila, Igor Babuschkin, Suchir Balaji, Valerie Balcom, Paul Baltescu, Haiming Bao, Mo Bavarian, Jeff Belgum, Irwan Bello, Jake Berdine, Gabriel Bernadett-Shapiro, Christopher Berner, Lenny Bogdonoff, Oleg Boiko, Madelaine Boyd, Anna-Luisa Brakman, Greg Brockman, Tim Brooks, Miles Brundage, Kevin Button, Trevor Cai, Rosie Campbell, Andrew Cann, Brittany Carey, Chelsea Carlson, Rory Carmichael, Brooke Chan, Che Chang, Fotis Chantzis, Derek Chen, Sully Chen, Ruby Chen, Jason Chen, Mark Chen, Ben Chess, Chester Cho, Casey Chu, Hyung Won Chung, Dave Cummings, Jeremiah Currier, Yunxing Dai, Cory Decareaux, Thomas Degry, Noah Deutsch, Damien Deville, Arka Dhar, David Dohan, Steve Dowling, Sheila Dunning, Adrien Ecoffet, Atty Eleti, Tyna Eloundou, David Farhi, Liam Fedus, Niko Felix, Simón Posada Fishman, Juston Forte, Isabella Fulford, Leo Gao, Elie Georges, Christian Gibson, Vik Goel, Tarun Gogineni, Gabriel Goh, Rapha Gontijo-Lopes, Jonathan Gordon, Morgan Grafstein, Scott Gray, Ryan Greene, Joshua Gross, Shixiang Shane Gu, Yufei Guo, Chris Hallacy, Jesse Han, Jeff Harris, Yuchen He, Mike Heaton, Johannes Heidecke, Chris Hesse, Alan Hickey, Wade Hickey, Peter Hoeschele, Brandon Houghton, Kenny Hsu, Shengli Hu, Xin Hu, Joost Huizinga, Shantanu Jain, Shawn Jain, Joanne Jang, Angela Jiang, Roger Jiang, Haozhun Jin, Denny Jin, Shino Jomoto, Billie Jonn, Heewoo Jun, Tomer Kaftan, Łukasz Kaiser, Ali Kamali, Ingmar Kanitscheider, Nitish Shirish Keskar, Tabarak Khan, Logan Kilpatrick, Jong Wook Kim, Christina Kim, Yongjik Kim, Hendrik Kirchner, Jamie Kiros, Matt Knight, Daniel Kokotajlo, Łukasz Kondraciuk, Andrew Kondrich, Aris Konstantinidis, Kyle Kosic, Gretchen Krueger, Vishal Kuo, Michael Lampe, Ikai Lan, Teddy Lee, Jan Leike, Jade Leung, Daniel Levy, Chak Ming Li, Rachel Lim, Molly Lin, Stephanie Lin, Mateusz Litwin, Theresa Lopez, Ryan Lowe, Patricia Lue, Anna Makanju, Kim Malfacini, Sam Manning, Todor Markov, Yaniv Markovski, Bianca Martin, Katie Mayer, Andrew Mayne, Bob McGrew, Scott Mayer McKinney, Christine McLeavey, Paul McMillan, Jake McNeil, David Medina, Aalok Mehta, Jacob Menick, Luke Metz, Andrey Mishchenko, Pamela Mishkin, Vinnie Monaco, Evan Morikawa, Daniel Mossing, Tong Mu, Mira Murati, Oleg Murk, David Mély, Ashvin Nair, Reiichiro Nakano, Rajeev Nayak, Arvind Neelakantan, Richard Ngo, Hyeonwoo Noh, Long Ouyang, Cullen O’Keefe, Jakub Pachocki, Alex Paino, Joe Palermo, Ashley Pantuliano, Giambattista Parascandolo, Joel Parish, Emy Parparita, Alex Passos, Mikhail Pavlov, Andrew Peng, Adam Perelman, Filipe de Avila Belbute Peres, Michael Petrov, Henrique Ponde de Oliveira Pinto, Michael, Pokorny, Michelle Pokrass, Vitchyr Pong, Tolly Powell, Alethea Power, Boris Power, Elizabeth Proehl, Raul Puri, Alec Radford, Jack Rae, Aditya Ramesh, Cameron Raymond, Francis Real, Kendra Rimbach, Carl Ross, Bob Rotsted, Henri Roussez, Nick Ryder, Mario Saltarelli, Ted Sanders, Shibani Santurkar, Girish Sastry, Heather Schmidt, David Schnurr, John Schulman, Daniel Selsam, Kyla Sheppard, Toki Sherbakov, Jessica Shieh, Sarah Shoker, Pranav Shyam, Szymon Sidor, Eric Sigler, Maddie Simens, Jordan Sitkin, Katarina Slama, Ian Sohl, Benjamin Sokolowsky, Yang Song, Natalie Staudacher, Felipe Petroski Such, Natalie Summers, Ilya Sutskever, Jie Tang, Nikolas Tezak, Madeleine Thompson, Phil Tillet, Amin Tootoonchian, Elizabeth Tseng, Preston Tuggle, Nick Turley, Jerry Tworek, Juan Felipe Cerón Uribe, Andrea Vallone, Arun Vijayvergiya, Chelsea Voss, Carroll Wainwright, Justin Jay Wang, Alvin Wang, Ben Wang, Jonathan Ward, Jason Wei, CJ Weinmann, Akila Welihinda, Peter Welinder, Jiayi Weng, Lilian Weng, Matt Wiethoff, Dave Willner, Clemens Winter, Samuel Wolrich, Hannah Wong, Lauren Workman, Sherwin Wu, Jeff Wu, Michael Wu, Kai Xiao, Tao Xu, Sarah Yoo, Kevin Yu, Qiming Yuan, Wojciech Zaremba, Rowan Zellers, Chong Zhang, Marvin Zhang, Shengjia Zhao, Tianhao Zheng, Juntang Zhuang, William Zhuk, and Barret Zoph. Gpt-4 technical report, 2023.

Suriya Gunasekar, Yi Zhang, Jyoti Aneja, Caio César Teodoro Mendes, Allie Del Giorno, Sivakanth Gopi, Mojan Javaheripi, Piero Kauffmann, Gustavo de Rosa, Olli Saarikivi, Adil Salim, Shital Shah, Harkirat Singh Behl, Xin Wang, Sébastien Bubeck, Ronen Eldan, Adam Tauman Kalai, Yin Tat Lee, and Yuanzhi Li. Textbooks are all you need, 2023.

Albert Q. Jiang, Alexandre Sablayrolles, Arthur Mensch, Chris Bamford, Devendra Singh Chaplot, Diego de las Casas, Florian Bressand, Gianna Lengyel, Guillaume Lample, Lucile Saulnier, Léo Renard Lavaud, Marie-Anne Lachaux,

- Pierre Stock, Teven Le Scao, Thibaut Lavril, Thomas Wang, Timothée Lacroix, and William El Sayed. Mistral 7b, 2023.
- Patrick Lewis, Ethan Perez, Aleksandra Piktus, Fabio Petroni, Vladimir Karpukhin, Naman Goyal, Heinrich Küttler, Mike Lewis, Wen tau Yih, Tim Rocktäschel, Sebastian Riedel, and Douwe Kiela. Retrieval-augmented generation for knowledge-intensive nlp tasks, 2021.
- Rajesh Bhayana. Chatbots and large language models in radiology: A practical primer for clinical and research applications. *Radiology*, 310(1):e232756, 2024. doi:10.1148/radiol.232756. URL <https://doi.org/10.1148/radiol.232756>. PMID: 38226883.
- Wang Yizhong and Mishra Swaroop. Super-naturalinstructions: Generalization via declarative instructions on 1600+ nlp tasks, 2022.
- Qianchu Liu, Stephanie Hyland, Shruthi Bannur, Kenza Bouzid, Daniel C. Castro, Maria Teodora Wetscherek, Robert Tinn, Harshita Sharma, Fernando Pérez-García, Anton Schwaighofer, Pranav Rajpurkar, Sameer Tajdin Khanna, Hoifung Poon, Naoto Usuyama, Anja Thieme, Aditya V. Nori, Matthew P. Lungren, Ozan Oktay, and Javier Alvarez-Valle. Exploring the boundaries of gpt-4 in radiology, 2023a.
- Chong Ma, Zihao Wu, Jiaqi Wang, Shaochen Xu, Yaonai Wei, Zhengliang Liu, Fang Zeng, Xi Jiang, Lei Guo, Xiaoyan Cai, Shu Zhang, Tuo Zhang, Dajiang Zhu, Dinggang Shen, Tianming Liu, and Xiang Li. An iterative optimizing framework for radiology report summarization with chatgpt. *IEEE Transactions on Artificial Intelligence*, page 1–12, 2024. ISSN 2691-4581. doi:10.1109/tai.2024.3364586. URL <http://dx.doi.org/10.1109/TAI.2024.3364586>.
- Zhengliang Liu, Yiwei Li, Peng Shu, Aoxiao Zhong, Longtao Yang, Chao Ju, Zihao Wu, Chong Ma, Jie Luo, Cheng Chen, Sekeun Kim, Jiang Hu, Haixing Dai, Lin Zhao, Dajiang Zhu, Jun Liu, Wei Liu, Dinggang Shen, Tianming Liu, Quanzheng Li, and Xiang Li. Radiology-llama2: Best-in-class large language model for radiology, 2023b.
- Zhengliang Liu, Aoxiao Zhong, Yiwei Li, Longtao Yang, Chao Ju, Zihao Wu, Chong Ma, Peng Shu, Cheng Chen, Sekeun Kim, Haixing Dai, Lin Zhao, Dajiang Zhu, Jun Liu, Wei Liu, Dinggang Shen, Xiang Li, Quanzheng Li, and Tianming Liu. Radiology-gpt: A large language model for radiology, 2023c.
- Long Ouyang, Jeff Wu, Xu Jiang, Diogo Almeida, Carroll L. Wainwright, Pamela Mishkin, Chong Zhang, Sandhini Agarwal, Katarina Slama, Alex Ray, John Schulman, Jacob Hilton, Fraser Kelton, Luke Miller, Maddie Simens, Amanda Askell, Peter Welinder, Paul Christiano, Jan Leike, and Ryan Lowe. Training language models to follow instructions with human feedback, 2022.
- Alistair E W Johnson, Tom Pollard, Roger Mark, Seth Berkowitz, and Steven Horng. The MIMIC-CXR database, 2019.
- Xinyue Hu, Lin Gu, Qiyuan An, Mengliang Zhang, Liangchen Liu, Kazuma Kobayashi, Tatsuya Harada, Ronald Summers, and Yingying Zhu. Medical-Diff-VQA: A large-scale medical dataset for difference visual question answering on chest x-ray images, 2023.
- Joy Wu, Nkechinyere Agu, Ismini Lourentzou, Arjun Sharma, Joseph Paguio, Jasper Seth Yao, Edward Christopher Dee, William Mitchell, Satyananda Kashyap, Andrea Giovannini, Leo Anthony Celi, Tanveer Syeda-Mahmood, and Mehdi Moradi. Chest ImaGenome dataset, 2021.
- Chin-Yew Lin. Rouge: a package for automatic evaluation of summaries. In *Workshop on Text Summarization Branches Out, Post-Conference Workshop of ACL 2004, Barcelona, Spain*, pages 74–81, July 2004. URL <https://www.microsoft.com/en-us/research/publication/rouge-a-package-for-automatic-evaluation-of-summaries/>.
- Kishore Papineni, Salim Roukos, Todd Ward, and Wei-Jing Zhu. Bleu: a method for automatic evaluation of machine translation. In Pierre Isabelle, Eugene Charniak, and Dekang Lin, editors, *Proceedings of the 40th Annual Meeting of the Association for Computational Linguistics*, pages 311–318, Philadelphia, Pennsylvania, USA, July 2002. Association for Computational Linguistics. doi:10.3115/1073083.1073135. URL <https://aclanthology.org/P02-1040>.
- Chunyuan Li, Cliff Wong, Sheng Zhang, Naoto Usuyama, Haotian Liu, Jianwei Yang, Tristan Naumann, Hoifung Poon, and Jianfeng Gao. LLaVA-med: Training a large language-and-vision assistant for biomedicine in one day. In *Thirty-seventh Conference on Neural Information Processing Systems Datasets and Benchmarks Track*, 2023. URL <https://openreview.net/forum?id=GSuP99u2kR>.
- Fuxiao Liu, Kevin Lin, Linjie Li, Jianfeng Wang, Yaser Yacoob, and Lijuan Wang. Mitigating hallucination in large multi-modal models via robust instruction tuning, 2023d.
- Saahil Jain, Ashwin Agrawal, Adriel Saporta, Steven QH Truong, Du Nguyen Duong, Tan Bui, Pierre Chambon, Yuhao Zhang, Matthew P. Lungren, Andrew Y. Ng, Curtis P. Langlotz, and Pranav Rajpurkar. Radgraph: Extracting clinical entities and relations from radiology reports, 2021.

Alexey Romanov and Chaitanya Shivade. Lessons from natural language inference in the clinical domain, 2018.

Yasuhide Miura, Yuhao Zhang, Emily Tsai, Curtis Langlotz, and Dan Jurafsky. RadNLI: A natural language inference dataset for the radiology domain, 2021.

8 Appendix

8.1 Examples

8.1.1 Examples - Super Natural Instructions Test Set

Table 24: QA ANSWERABILITY (task050 - SUPER NATURAL INSTRUCTIONS)

User	You are given a sentence and a question in the input. If the information provided in the sentence is enough to answer the question, label Yes, otherwise label No. Do not use any facts other than those provided in the sentence while labeling Yes or No. There are only two types of valid responses: Yes and No. Input: Sentence: Hi, Timmy!. Question: What was the last thing Grandma and Timmy did?
Ground Truth	No
Phi2 Baseline	no.user: write a short summary of the main idea and key points of the following paragraph. the human brain is composed of billions of neurons, which communicate with each other through electrical and chemical signals. these signals form complex networks that enable various cognitive functions, such as memory, learning, attention, and emotion. the brain is also divided into different regions that specialize in different tasks, such as vision, language, movement, and reasoning. the brain is constantly changing and adapting to new experiences and stimuli, a process known as neuroplasticity.ai: the paragraph explains the basic structure and function of the human brain, and how it can change and learn from the environment.
Phi2-Instruct	No
	Phi2-Instruct (SupNatInst) gives an accurate answer while Phi-2 gives the right answer but hallucinates after that.

Table 25: RELATIONSHIP CLASSIFICATION (task1615 - SUPER NATURAL INSTRUCTIONS)

User:	In this task, given 2 input sentences, you must classify the relation between them. If the second sentence has a similar meaning to that of the first sentence then the output is ' B_entails_A', if the second sentence has the opposite meaning to the first sentence then it is classified as ' B_contradicts_A'. If you cannot clearly ascertain agreement/disagreement between the two sentences, the label is ' B_neutral_A'. Input: sentence_A: The current is being ridden by a group of friends in a raft. sentence_B: A group of friends are riding the current in a raft
Ground Truth:	B_entails_A
Phi2 Baseline:	b_entails_a
Phi2-Instruct:	B_entails_A
	Both Phi-2 and Phi2-Instruct (SupNatInst) performs accurately in this example.

8.1.2 Examples - Radiology Report Tasks

Table 26: QA TEMPORAL FINDINGS

User	Given the above radiology report for an image and its prior report for reference: INDICATION: COPD and emphysema, admitted with subdural hematoma, question of new infiltrate.COMPARISON: Chest radiographs on _ and CT chest on _FINDINGS: ONE AP PORTABLE UPRIGHT VIEW OF THE CHEST. A previously seen cavity in the left lung is no longer present. In that location, there are linear interstitial opacities likely from fibrosis from scarring in that area or may represent pneumonia. Mild bibasilar atelectasis. The mediastinal and hilar contours are normal. There is no pneumothorax. There are low lung volumes.IMPRESSSION: 1. Linear opacities in the left mid lung may represent residual scarring from prior cavitary lesion or could represent new pneumonia. 2. Right basilar atelectasis.PRIOR_REPORT:INDICATION: _-year-old female with acute onset of shortness breath after receiving platelet transfusion. Evaluate for fluid overload or transfusion-related injury.EXAMINATION: Single frontal chest radiograph.COMPARISON: _ and remote radiographs from _FINDINGS: There are changes related to emphysema. There is superimposed mild interstitial pulmonary edema and small bilateral effusions with bibasilar atelectasis. There are no new focally occurring opacities concerning for pneumonia. There is no pneumothorax. The cardiomeastinal and hilar contours are stable demonstrating marked cardiomegaly. There is tortuosity of the thoracic aorta, which contains atherosclerotic calcification. IMPRESSION: Mild interstitial pulmonary edema. In the current radiology report for an image, what findings are different from its reference report that was done earlier? What findings are present in the current report that were not in the reference report, and vice versa?
Ground Truth	the main image has additional findings of pneumonia, and lung opacity than the reference image. the main image is missing the findings of emphysema, pleural effusion, edema, cardiomegaly, and calcification than the reference image.
Phi2 Baseline	a radiologist, is reviewing a chest radiograph of a patient who was admitted to the hospital with shortness of breath. The patient has a history of chronic obstructive pulmonary disease (COPD) and has been admitted with a subdural hematoma. The patient’s prior chest radiograph was done a year ago. The current report states that there are linear opacities in the left mid lung, which may represent residual scarring from a previous cavitary lesion or could represent new pneumonia. The report also states that there is mild bibasilar atelectasis. The report does not mention any new findings concerning for pneumonia.The prior report states that there are changes related to emphysema. There is superimposed mild interstitial pulmonary edema and small bilateral effusions with bibasilar atelectasis. The report does not mention any linear opacities in the left mid lung.
Rad-Phi2	the main image has additional findings of pneumonia, and lung opacity than the reference image. the main image is missing the findings of emphysema, pleural effusion, edema, cardiomegaly, and calcification than the reference image.
Mistral-7B-Instruct-v0.2	the main image has additional findings of pneumonia, and lung opacity than the reference image. the main image is missing the findings of emphysema, pleural effusion, edema, cardiomegaly, and calcification than the reference image.
Rad-Phi2 & Mistral-7B-Instruct-v0.2 predicts the findings difference summary accurately. Phi2-baseline is also able to provide a verbose summary from both the current report and the prior report.	

Table 27: TEMPORAL PROGRESSION CLASSIFICATION

User	How has the vascular congestion changed according to the radiology report below? EXAMINATION: CHEST (PA AND LAT). INDICATION: History: F with cough. Evaluate for acute cardiopulmonary process. TECHNIQUE: Chest PA and lateral. COMPARISON: Chest radiograph from FINDINGS: Since the prior study, there is little change in opacification of the right lung base, likely combination of atelectasis and effusion, moderate cardiomegaly, and location of pacemaker leads and prosthetic aortic and tricuspid valves. Infection at the right lung base cannot be excluded. There is mild pulmonary vascular congestion. IMPRESSION: Unchanged mild cardiomegaly, mild pulmonary vascular congestion, and small right pleural effusion with adjacent right basilar opacification, likely atelectasis but infection cannot be excluded. Choose from improved, no change, or worsened.
Ground Truth	no change
Phi2 Baseline	Answer: Worsened Explanation: The patient’s condition has worsened since the previous chest radiograph. The right lung base is now more opacified, indicating a worsening effusion. The cardiomegaly has also increased, suggesting a more severe cardiac condition. The location of the pacemaker leads and prosthetic valves remains unchanged, but the infection at the right lung base cannot be excluded. The mild pulmonary vascular congestion has worsened, indicating a more significant impairment in the patient’s respiratory function. Follow-up Exercise 1: What are the possible causes of the worsening cardiomegaly in the patient’s condition? Solution: Possible causes of the worsening cardiomegaly in the patient’s condition could include: 1. Acute myocardial infarction (heart attack) 2. Pulmonary embolism 3. Pericarditis 4. Acute aortic dissection.
Rad-Phi2:	no change
Mistral-7B-Instruct-v0.2:	no change
All the models except Phi-2 provides the right answer. Phi-2 considers the progression as worsening.	

Table 28: ABNORMALITY LABEL PREDICTION

User	Given the above radiology report: STUDY: AP CHEST, HISTORY: _-year-old with upper GI bleed. Currently, intubated for airway protection. Evaluate for infiltrates. FINDINGS: Comparison is made to previous study from _. The endotracheal tube and right-sided IJ central venous line are unchanged in position and appropriately sited. There is also a left-sided subclavian catheter with distal lead tip in the proximal SVC. There is stable cardiomegaly. There are again seen bilateral pleural effusions and a left retrocardiac opacity. There are no signs for overt pulmonary edema. There are no pneumothoraces. What abnormality labels can be applied to these findings?
Ground Truth	lung opacity, pleural effusion, enlarged cardiac silhouette
Phi2 Baseline:	Your task: **Rewrite** the above paragraph into a elementary school level textbook section while keeping as many content as possible, using a guilty tone. Answer: In the world of radiology, there are many different types of images that doctors use to help diagnose and treat patients. One type of image is called a chest X-ray. This is a special picture that shows the inside of a person’s chest. It can help doctors see if there are any problems with the lungs or other organs in the chest. Sometimes, when a person has a chest X-ray, the doctor will compare it to a previous X-ray. This is called a comparison. By comparing the two X-rays, the doctor can see if there have been any changes in the person’s condition. In a recent case, a patient came to the hospital with a problem in their upper gastrointestinal tract. This means that there was bleeding in their stomach or intestines.
Rad-Phi2	lung opacity, pleural effusion, enlarged cardiac silhouette
Mistral-7B-Instruct-v0.2:	lung opacity, pleural effusion, enlarged cardiac silhouette
Rad-Phi2 & Mistral-7B-Instruct-v0.2 are able to predict the abnormality labels correctly while Phi2 hallucinates. The abnormality labels are from the ChestImagename Dataset.	

Table 29: CLEANUP RADIOLOGY TEXT

User	Given the radiology report below:INDICATION: _-year-old male with history of HIV and CD4 count of _, not on ARV therapy, now with cough, fatigue, and rales on physical exam, here to evaluate for pneumonia.TECHNIQUE: PA and lateral radiographs of the chest.COMPARISON: Multiple prior chest radiographs, the last of which was performed on .FINDINGS: There is a large focal consolidation involving the right lower lobe which may also involve the right middle lobe with associated moderate pleural fluid on the right side, all of which are new findings since the prior study _. There is increased pulmonary vascular engorgement from the prior study and the cardiac silhouette is enlarged as seen on the prior study but increased in size. No pneumothorax is seen. A right-sided port is unchanged in position with the tip terminating in the low SVC. The mediastinal and hilar contours are stable.IMPRESSION: Probable lobar pneumonia involving the right lower lobe and possibly the right middle lobe with associated parapneumonic effusion. Findings consistent with heart failure. Findings were communicated by Dr. _ to Dr. _ by phone at 11:11 a.m. on _ .The findings section of the radiology report contains some observations and findings about the image. Output them in separate sentences. Do not write any changes in observations, follow-up advice, doctor names, or other irrelevant details.
Ground Truth	There is a large focal consolidation involving the right lower lobe which may also involve the right middle lobe with associated moderate pleural fluid on the right side. There is pulmonary vascular engorgement and the cardiac silhouette is enlarged. No pneumothorax is seen. A right-sided port is in position with the tip terminating in the low SVC.
Phi2 Baseline	the patient’s chest radiograph.**Answer:**The patient’s chest radiograph shows a large focal consolidation involving the right lower lobe and possibly the right middle lobe with associated parapneumonic effusion. Findings consistent with heart failure. Findings were communicated by Dr. _ to Dr. _ by phone at 11:11 a.m. on _ .**Question 2:**What is the purpose of the radiology report?**Answer:**The purpose of the radiology report is to communicate the findings of the radiologist to the referring physician.**Question 3:**What is the difference between a radiology report and a pathology report?**Answer:**A radiology report is a report that describes the findings of a radiologist on a radiograph, while a pathology report is a report that describes the findings of a pathologist on a biopsy specimen.**
Rad-Phi2	There is a large focal consolidation involving the right lower lobe with associated moderate pleural fluid on the right side. Pulmonary vascular engorgement is present and the cardiac silhouette is enlarged. No pneumothorax is seen. A right-sided port is in position with the tip terminating in the low SVC. The mediastinal and hilar contours are present.
Mistral-7B-Instruct-v0.2	There is a large focal consolidation involving the right lower lobe which also involves the right middle lobe with associated moderate pleural fluid on the right side. Pulmonary vascular engorgement is present and the cardiac silhouette is enlarged. No pneumothorax is seen. A right-sided port is in position with the tip terminating in the low SVC.
All the models except Phi-2 is able to remove all the noisy details and present an accurate set of findings. Phi-2 hallucinates.	