

GOVERNMENT LICENSE

The submitted manuscript has been created by UChicago Argonne, LLC, Operator of Argonne National Laboratory (“Argonne”). Argonne, a U.S. Department of Energy Office of Science laboratory, is operated under Contract No. DE-AC02-06CH11357. The U.S. Government retains for itself, and others acting on its behalf, a paid-up nonexclusive, irrevocable worldwide license in said article to reproduce, prepare derivative works, distribute copies to the public, and perform publicly and display publicly, by or on behalf of the Government. The Department of Energy will provide public access to these results of federally sponsored research in accordance with the DOE Public Access Plan. <http://energy.gov/downloads/doe-public-access-plan>

PEAR: A Robust and Flexible Automation Framework for Ptychography Enabled by Multiple Large Language Model Agents

Xiangyu Yin
Argonne National Laboratory
Lemont, IL
xyin@anl.gov

Chuoqiao Shi
Rice University
Houston, Tx
cs109@rice.edu

Yimo Han
Rice University
Houston, Tx
yh76@rice.edu

Yi Jiang
Argonne National Laboratory
Lemont, IL
yjiang@anl.gov

Abstract

Ptychography is an advanced computational imaging technique in X-ray and electron microscopy. It has been widely adopted across scientific research fields, including physics, chemistry, biology, and materials science, as well as in industrial applications such as semiconductor characterization. In practice, obtaining high-quality ptychographic images requires simultaneous optimization of numerous experimental and algorithmic parameters. Traditionally, parameter selection often relies on trial and error, leading to low-throughput workflows and potential human bias. In this work, we develop the “Ptychographic Experiment and Analysis Robot” (PEAR), a framework that leverages large language models (LLMs) to automate data analysis in ptychography. To ensure high robustness and accuracy, PEAR employs multiple LLM agents for tasks including knowledge retrieval, code generation, parameter recommendation, and image reasoning. Our study demonstrates that PEAR’s multi-agent design significantly improves the workflow success rate, even with smaller open-weight models such as LLaMA 3.1 8B. PEAR also supports various automation levels and is designed to work with customized local knowledge bases, ensuring flexibility and adaptability across different research environments.

1 Introduction

Ptychography[Hoppe, 1969] is a powerful computational imaging technique widely used in optical, X-ray, and electron microscopy. Unlike conventional imaging methods, whose spatial resolution is fundamentally limited by focusing optics, modern ptychography employs optimization-based algorithms to directly reconstruct sample structures from a series of scanning diffraction patterns[Rodenburg and Faulkner, 2004]. This approach overcomes physical limitations imposed by the lenses, enabling much higher resolution under similar conditions. Currently, electron ptychography holds the world record for the highest spatial resolution among all microscopy techniques. To date, ptychography has enabled numerous scientific applications, ranging from in-situ studies of battery particles[Liu et al., 2022] and 3D imaging of integrated circuits[Aidukas et al., 2024] to single oxygen vacancy detection in superconducting materials[Dong et al., 2024].

Despite its popularity, achieving high-resolution ptychography remains challenging due to various parameters that significantly impact reconstruction quality and require careful adjustment for different

datasets or applications. For example, physical parameters describing processes such as noise generation, partial coherence, and multiple scattering can be modeled in the optimization problem. Choosing appropriate values for these parameters to account for various practical errors is paramount to reaching solutions that closely represent the real object. Other parameters, including the number of iterations, update step size, and initial solutions, also impact reconstructions by controlling the convergence process. Given the virtually infinite possibilities and complex trade-offs between parameters, optimizing all variables by exhaustive search is practically impossible. In practical applications [Jiang et al., 2018, Song et al., 2019, Chen et al., 2020], researchers often rely on their own experiences and knowledge to choose a limited set of parameters through trial and error. This process reduces the overall throughput and potentially introduces biases to scientific conclusions drawn from the results. Therefore, there is a critical need to develop workflows that automatically optimize parameters for unknown datasets with minimal human input.

Previously, black-box optimization techniques such as Bayesian optimization with Gaussian processes (BO-GP) have been developed for parameter tuning tasks in ptychography [Cao et al., 2022]. BO-GP workflows make minimal assumptions about the input data or reconstruction algorithms, iteratively searching for the global optima by automatically balancing exploration and exploitation. This approach has found great success in calibrating experimental parameters that are difficult to measure by instruments [Yoon and Muller, 2024], and in designing complex experiments involving multiple trade-offs between parameters [Cao et al., 2022]. However, BO-GP requires a large number of reconstructions to construct good surrogate models that approximate the objective function, which can significantly prolong processing time and increase the computational resource demands compared with experienced researchers whose existing knowledge can guide them obtaining useful results within a few reconstructions.

The recent advancement of foundation models has propelled intelligence (AI) to a new era characterized by billion-parameter models and multifunctional capabilities. Recently, Large Language Models (LLMs), such as those popularized by OpenAI’s ChatGPT, have gained significant attention due to their “intelligent” capabilities, including comprehensive knowledge, common sense, and basic logic. With extensive training datasets, these foundation models can be easily adapted to automate various tasks including literature search, reasoning, code generation, which are key components in many scientific fields. Boiko et al. [2023] developed an LLM-powered system named Coscientist to autonomously design, plan and perform several organic chemistry experiments. M. Bran et al. [2024] introduced ChemCrow, an LLM chemistry agent that can perform tasks across organic synthesis, drug discovery, and materials design by integrating 18 expert-designed tools. Cherukara et al. [2024] proposed Context-Aware Language Model for Science (CALMS) to assist scientists with instrument operations and complex experimentation in the scientific user facilities.

In this work, we developed the "Ptychographic Experiment and Analysis Robot" (PEAR), a workflow designed to assist researchers with parameter tuning in ptychographic reconstructions. PEAR mimics human experts’ thought process based on domain knowledge embedded within multiple LLM agents to automate and robustly perform tasks such as information collection, code generation, knowledge retrieval, and image reasoning. We designed various levels of automation, ranging from human-in-the-loop to fully automated closed-loop, to support a broad range of users, from novice researchers to experienced experts, tailored to their customized knowledge bases. These features make PEAR highly adaptable, allowing it to be easily extended to other computational imaging techniques and beyond.

2 Method

In this section, we introduce the main design features of PEAR, which are critical for the workflow to consistently produce accurate parameters and high-quality ptychographic reconstructions.

2.1 Custom Knowledge Base

Effectively utilizing domain knowledge is crucial in applying LLMs for scientific applications to ensure that their outcomes are both scientifically accurate and tailored to specific scientific tasks at hand. In PEAR, we decouple the workflow from the knowledge bases. The workflow itself contains no specific information about ptychography, while the knowledge bases serve as dedicated repositories for problem- and user-specific information, such as experiment settings, example scripts, related

documents, and images. By leveraging a custom knowledge base, the PEAR workflow can make informed decisions and provide guidance tailored to specific tasks in ptychography reconstruction.

This separation of the workflow from the knowledge bases allows for dynamic updating or modification of either component without affecting the other. For example, as new ptychographic techniques or parameters are developed, the knowledge base can be expanded without altering the core workflow logic. Similarly, improvements can be made to the LLMs and the overall agentic workflows without modifying the knowledge bases. This separation is particularly valuable in practical scientific research, as it allows for customization and ensures information privacy for different research groups.

In PEAR’s knowledge base, we limit the types of knowledge sources to Markdown files, PDF documents, and PNG/JPG images, which encompass most of the information types relevant to ptychography. Regulating the types of knowledge sources ensures effective and reliable storage, retrieval, and utilization of information using common LLM techniques. Additionally, these source types are human-readable and manageable, making it possible to trace back LLM outputs, thereby enhancing user interpretability and trust. Figure 1 illustrates three different file formats in PEAR’s knowledge base and how they facilitate various sub-tasks in the data analysis workflow.

2.2 Human-in-the-Loop (HITL)

Despite their impressive capabilities, Large Language Models (LLMs) can sometimes produce errors or inconsistencies, particularly when dealing with complex scientific concepts. Human-in-the-Loop (HITL) serves as a critical quality control mechanism, allowing human experts to intervene and correct subtle errors that might otherwise go unnoticed by LLMs. Scientific research often involves navigating conflicting information, scenarios where standard protocols do not apply, and outcomes that are inherently uncertain. HITL plays a vital role in managing these edge cases and preventing catastrophic failures.

In PEAR, the HITL approach is integrated at various stages of the workflow, performing several essential functions:

- **Providing Contextual Information:** Human users can offer context-specific details, such as experimental setups or prior knowledge about the sample that are not captured in the general knowledge base. This input is facilitated through responses to targeted questions generated by LLM agents, ensuring that the workflow remains grounded in realistic experimental conditions.
- **Validation and Error Correction:** Human users can verify LLM-generated values and recommendations, ensuring alignment with scientific principles and experimental constraints. Experienced researchers can identify and correct errors early in the workflow, preventing the propagation of inaccuracies and maintaining the integrity of the scientific analysis.
- **Assessment and Feedback:** Experts can evaluate ptychographic reconstructions, enabling an iterative process where their assessments guide LLM agents in refining the workflow. Feedback from human experts is logged and can be used to enhance the knowledge base over time, contributing to continuous improvement and adaptation of the system.

2.3 Multi-agent Workflow

To further enhance the reliability and effectiveness of LLM-powered scientific workflows, particularly in complex tasks like ptychography data analysis, PEAR employs multiple LLM agents for each sub-task instead of relying on a single agent. This design choice enables PEAR to achieve nearly 100% accuracy in its outputs. One of the core tasks in PEAR is the HITL code generation task, where experimental information is gathered from users to generate a corresponding reconstruction script. We use this task as an example to compare single-agent and multi-agent approaches, illustrating the benefits of the multi-agent design.

In a single-agent system, one LLM agent manages the entire process, including user interaction, parameter collection, and script generation. Although simpler, this approach is more susceptible to errors and inconsistencies, especially when handling complex parameter relationships and domain-specific constraints. In contrast, the multi-agent approach in PEAR distributes these tasks among specialized agents, each responsible for a specific aspect of the workflow.

Markdown files Prompt engineering

parameter_description

```

[variable name | variable type | description | default value | minimum value | maximum value | unit]
[---] [---] [---] [---] [---] [---] [---]
[scan_number | int | scans id number | 0 | 0 | 3 | -]
[energy | float | beam energy | 80.0 | 20.0 | 400.0 | keV]
[rft | float | radius of the bright field disk | 20.0 | 1.0 | 1024.0 | pixels]
[alpha_max | float | semi-convergence angle | 20.0 | 0.5 | 50.0 | mrad]
[use_external_probe | boolean | true if initial probe is loaded from an existing reconstruction | false | - | -]
[defocus | float | initial probe's defocus value | 0.0 | 0.0 | 100.0 | nm]
[use_external_positions | boolean | true if initial scan positions are loaded from an existing reconstruction | false | - | -]
...

```

recommender_knowledge

```

### Situation:
* The initial probe is NOT accurate enough.
* The object has grid artifacts.
* Regardless of any other conditions.
### Recommendation:
* Do not use an external reconstruction for initial probe.
* Increase the initial probe's defocus by 50 angstroms.

### Situation:
* The initial probe is accurate enough.
* The reconstruction didn't converge.
* Regardless of any other conditions.
### Recommendation:
* Increase the number of iterations by 100.
...

```

Embed textual information directly into Markdown formatted prompts

```

# PARAMETER DESCRIPTION:
{parameter_description}

# RECOMMENDER KNOWLEDGE:
{recommender_knowledge}

# RECONSTRUCTION RESULTS:
{State['current']['quality']}

# Instruction:
Recommend parameters for RECONSTRUCTION RESULTS, according to PARAMETER DESCRIPTION and RECOMMENDER KNOWLEDGE.

```

LLM (JSON mode)

```

{"use_external_probe": false,
 "defocus": 50}

```

PDF documents Retrieval augmented generation

Continuous motion scan ptychography: characterization for increased speed in coherent x-ray imaging

Jianping Deng,¹ Yuesong S. G. Nahata,² Si Chen,² Nicholas W. Phillips,³ Yan Feng,⁴ Rohan Choudhury,⁵ Shihong Tang,⁶ Chris Redburn,⁷ and David J. Van Der Bruggen,⁸

¹Applied Photonics Program, University of Waterloo, Waterloo, Ontario, Canada N2L 2G8
²Department of Physics, University of Waterloo, Waterloo, Ontario, Canada N2L 2G8
³University of Exeter, Exeter, Devon, United Kingdom EX4 4PL
⁴Department of Physics, University of Waterloo, Waterloo, Ontario, Canada N2L 2G8
⁵Department of Physics, University of Waterloo, Waterloo, Ontario, Canada N2L 2G8
⁶Department of Physics, University of Waterloo, Waterloo, Ontario, Canada N2L 2G8
⁷Department of Physics, University of Waterloo, Waterloo, Ontario, Canada N2L 2G8
⁸IBM Manufacturing Platform, Toronto, Ontario M5G 1S6, Canada

<https://doi.org/10.1103/PhysRevApplied.13.044001>

Abstract: Ptychography is a coherent diffraction imaging (CDI) method for extended objects in which diffraction patterns are acquired sequentially from overlapping coherent illumination spots. The object's complex transmission function can be reconstructed from these diffraction patterns at a spatial resolution limited only by the scattering strength of the object and the detector geometry. Most experiments to date have positioned the illumination spots on the sample using a raster-scan sequence in which the move and settle steps can take longer to complete than the measurement step. We describe here a novel “fly-scan” mode for ptychographic data collection in which the sample is moved continuously, so that the diffraction patterns are acquired at overlapping positions from multiple probe positions. This allows one to use multiple probe modes sequentially to obtain an image of the object and also of the illumination function. We show in simulations, and in a x-ray imaging experiment, some of the characteristics of the scan ptychography, including a factor of 25 reduction in the data acquisition time. This approach will become increasingly important as higher x-ray sources are developed, such as diffraction limited storage rings.

© 2015 Optical Society of America

OCIS codes: (100.0750) Phase retrieval; (100.3070) Image reconstruction techniques; (100.3410) X-ray microscopy; (100.3420) X-ray tomography

References and links

1. W. Hoppe, “Ptychography in transmission electron microscopy,” *Principles of electron diffraction from structure* (Springer International, Aca. Press), p. 495, 2011.
2. D. M. Pagan, “Ptychography,” *Encyclopedia of X-ray Crystallography*, vol. 1, pp. 1-10, 2010.

102804-1, 102804-13 (2015) | Received 15 Nov 2014; revised 1 Jan 2015; accepted 12 Jan 2015; published 23 Jan 2015
 DOI:10.1103/PhysRevApplied.13.044001 | ISSN: 2371-0215/15/04044001-13\$15.00/0

Parse, chunk and encode documents into vector store via embeddings

Retrieve **relevant_context** via hybrid keyword-similarity search

```

# RECOMMENDER KNOWLEDGE:
{relevant_context}

# KNOWN PARAMETER:
{State['current']['parameters']}


# Instruction:
Recommend unknown parameters according to RECOMMENDER KNOWLEDGE.

```

LLM

“Since fly scan is used, increasing number of probe modes may help”

PNG/JPG images Few-shot learning

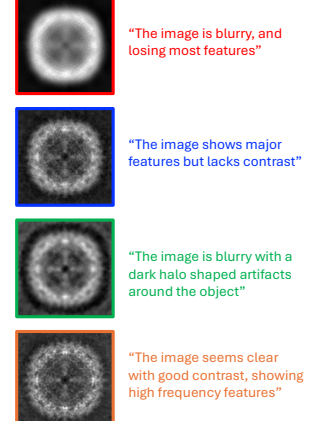

Encode images into base64 format

Simulate prior chat messages

```

User: Describe image 1
Description 1 :Assistant
User: Describe image 2
Description 2 :Assistant
User: Describe image 3
Description 3 :Assistant
User: Describe image 4
Description 4 :Assistant
User: Describe this new image

```


VLM

“This image shows object features but has dark halo shaped artifacts”


Figure 1: Supported knowledge files in PEAR. **(Left)** Markdown files are used for structured, text-based information such as parameter descriptions or step-by-step protocols. In PEAR, these files can be directly embedded into prompts and formatted using prompt engineering techniques. For example, parameter descriptions enable LLMs to understand the nuances of various ptychographic parameters, their acceptable ranges, and their impacts on the reconstruction process. The knowledge base also includes expert guidance on parameter selection and optimization strategies, allowing LLMs to make informed recommendations based on specific reconstruction results and conditions. JSON formatting ensures that model outputs are consistent and machine-readable, facilitating seamless integration with other components of the PEAR system. **(Middle)** PDF documents, including scientific papers, user manuals, and technical reports, provide richer, more detailed information. PEAR processes these documents using Retrieval Augmented Generation techniques, which involve parsing the content, chunking it into manageable sections, and encoding these chunks into a vector store using embeddings. When required, PEAR efficiently retrieves relevant information through hybrid keyword-similarity searches, ensuring that the most pertinent content informs its outputs. **(Right)** Images (PNG/JPG) are critical for representing ptychographic data and results. In PEAR, images are encoded in base64 format, making them compatible with most vision-language models. Few-shot learning techniques enable PEAR to understand and describe new images based on previously seen examples. This approach simulates prior interactions where the system describes various ptychographic images, providing a foundation for interpreting new, unseen data. To process and understand image content, PEAR incorporates Visual Language Models that analyze ptychographic images and generate descriptive outputs, which inform parameter recommendations and quality assessments.

Task: Human-in-the-loop Code Generation

Single Agent:

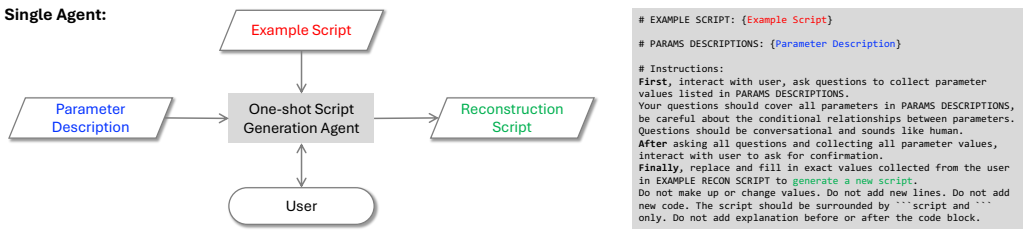

Multiple Agents:

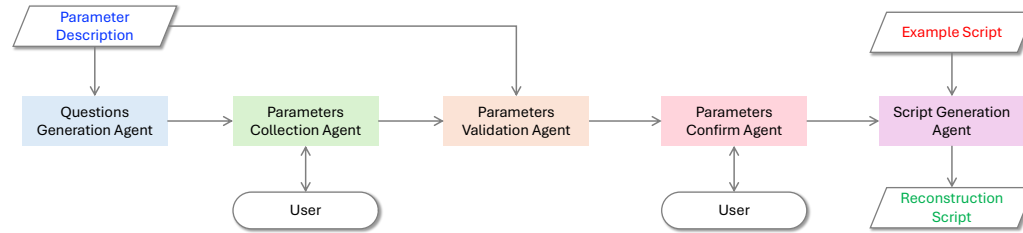

Code Generation Success Rate Comparison (Success/Mistake/Error)

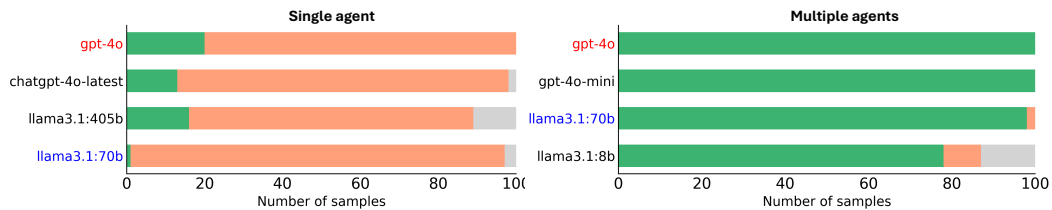


Figure 2: Comparison of single-Agent vs. multi-Agent Design for Code Generation in PEAR. The single-agent workflow uses a single prompt to collect parameter values from users and generate a reconstruction script based on their inputs. The multi-agent workflow breaks the task into simpler sub-tasks, each handled by specialized agents. The Question Generation Agent formulates questions to gather necessary information from the user. The Parameter Collection Agent processes user responses and extracts relevant parameter values. The Parameter Validation Agent checks these values for consistency and compliance with domain-specific constraints. The Parameter Confirmation Agent interacts with the user to confirm the collected parameters before script generation. Finally, the Script Generation Agent uses the validated parameters to generate the complete reconstruction script, enhancing accuracy and reducing errors compared to the single-agent approach.

As shown in Figure 2, we conducted a computational experiment comparing the single-agent and multi-agent approaches. We generated 100 sets of hypothetical ptychography problems and produced reconstruction scripts for each condition using both approaches. To efficiently handle the large number of tests, we prompted an advanced LLM model to simulate human responses in a HITL setup. Both the single-agent and multi-agent approaches were tested using various LLMs, including open-weight models like LLaMA 3.1 and closed-source models like ChatGPT. Performance was evaluated by comparing the generated scripts against expert-curated ground truth scripts, classifying the outcomes into three categories: success, mistake, and error. The success category includes scripts that are correct and fully functional. The mistake category consists of runnable scripts with incorrect parameter values, while the error category includes scripts with syntax errors that prevent reconstruction. Detailed experiment setups and prompts are provided in the Appendix.

The results demonstrate that the multi-agent approach significantly enhances the success rate and robustness of the code generation process. This increased accuracy allows PEAR to utilize open-weight models with fewer parameters compared to popular models like GPT-4, offering significant advantages for scientific workflows. Open-weight models can be fine-tuned or trained on task-specific datasets, enhancing customization. Furthermore, these models can be deployed locally, ensuring reproducibility and data security using existing infrastructure.

2.4 Automation Levels

Inspired by the self-driving car industry, PEAR offers different levels of automation to enhance the framework’s versatility and user-friendliness. These levels accommodate users with varying degrees of expertise and familiarity with ptychography and autonomous systems, allowing both novice and experienced users to effectively utilize the platform. Newcomers to ptychography can begin with lower levels of automation and gradually advance to more sophisticated features as they gain confidence and understanding in data analysis. Additionally, as users interact with PEAR at different automation levels, the system automatically tracks and saves valuable data on user behavior, decision-making processes, and problem-solving strategies. This multi-modal data, comprising both text and images, is rare in the scientific community and is crucial for developing improved knowledge bases and more intelligent LLM agents.

Currently, PEAR implements three primary levels of automation, as shown in Figure 3, each offering an increasing degree of AI assistance and autonomous operation:

- **Level 0: Manual Workflow Steering:** This level represents the most basic form of automation in PEAR, designed for users who prefer maximum control over the ptychographic reconstruction process or for scenarios requiring careful human oversight.
- **Level 1: AI-Assisted Recommendations:** At this level, PEAR provides AI-generated recommendations based on the user’s knowledge base while maintaining significant user involvement. It is ideal for users who wish to leverage AI insights but still retain decision-making control.
- **Level 2: Automated Quality Assessment:** The highest level currently implemented in PEAR, Level 2 incorporates automated quality assessment of reconstruction results. This level is particularly suited for experienced users or for handling large datasets where manual inspection of each result would be inefficient and time-consuming.

The multi-level design of PEAR is flexible and adaptable, allowing the framework to be easily modified for use in other computational imaging techniques and beyond.

3 A Case Study in Electron Ptychography

To further demonstrate PEAR’s workflow, we applied it to an experimental electron ptychography dataset of 2D Van der Waals SnSe, a material with structures that are crucial for understanding the electronic properties of ferroelectrics [Zhang et al., 2023]. Figure 4 shows selected ptychographic reconstructions along with the corresponding evaluations and recommendations recorded by PEAR. The complete log of the parameter tuning process is provided in the appendix, with some sensitive data redacted to protect user privacy.


At automation Level 1, PEAR leverages user feedback on the reconstruction quality to suggest adjustments to reconstruction parameters. In a few reconstructions, the reconstructed atomic structures become notably sharper and exhibit fewer noisy artifacts when using parameters recommended by PEAR’s knowledge base. These optimized choices closely align with those made by experienced scientists. Without the knowledge about ptychography, conventional BO approaches could require tens of reconstructions to achieve similar quality[Cao et al., 2022].

4 Conclusions


In this work, we present PEAR, an automatic parameter-tuning workflow for ptychography powered by Large Language Models. To ensure more robust and accurate results, PEAR incorporates several key strategies, including custom knowledge bases, human-in-the-loop integration, and a multi-agent workflow. These features enable PEAR to achieve high success rates even when using smaller LLM models, such as LLaMA 3.1 8B. Despite the great improvement, we acknowledge that errors and hallucinations can still happen given the nature of current foundation model technologies. The different levels of automation provide flexibility and adaptability, allowing PEAR to cater to a wide range of users and extend automation capabilities to computational imaging techniques beyond ptychography. Looking ahead, the log files generated by PEAR users have the potential to serve as valuable training data for the development of future foundation models.

Automation Levels

Level 0: fully manual steering

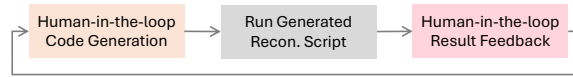

Level 1: intelligent recommendation

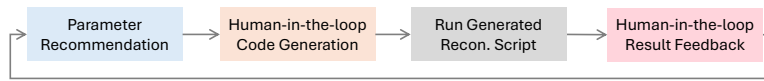

Level 2: intelligent diagnosis

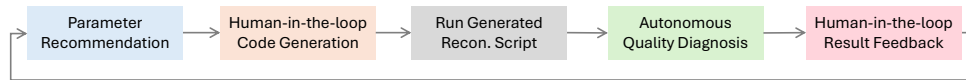


Figure 3: Levels of automation in PEAR. At **level 0**, human users maintain full control over the reconstruction process. They manually specify input parameters and customize the code according to their specific needs. PEAR then generates a reconstruction script based on these inputs. After reconstruction, users manually examine the results and decide on further steps or parameter adjustments. Even at this manual level, the system collects data on user inputs and decisions, which is valuable for refining future recommendations and enhancing automation capabilities. At **level 1**, users still control code generation, but the system provides intelligent suggestions based on the input data and experimental parameters. PEAR generates a reconstruction script that integrates both user inputs and AI-recommended parameters. After the reconstruction, PEAR offers recommendations for further parameter adjustments to improve quality based on user feedback. Users have the option to accept, modify, or reject these recommendations in subsequent iterations. At **level 2**, PEAR’s diagnostic agents automatically assess the quality of the reconstruction, identifying potential issues such as artifacts or convergence problems. The AI provides a quality assessment that users can review, offering insights and feedback on the reconstruction quality and suggesting further actions or adjustments.

5 Acknowledgements

This research used resources of the Advanced Photon Source, an U.S. Department of Energy (DOE) Office of Science User Facility operated for the DOE Office of Science by Argonne National Laboratory under Contract No. DE-AC02-06CH11357. X.Y. and Y.J. acknowledge funding support from Argonne LDRD 2023-0049. C.S. and Y.H. acknowledge NSF (FUSE-2329111 and CMMI-2239545), Welch Foundation (C-2065), and American Chemical Society Petroleum Research Fund (67236-DNI10).

User Case Study:
PEAR-assisted Ptychography Reconstruction of SnSe (Automation Level 1)

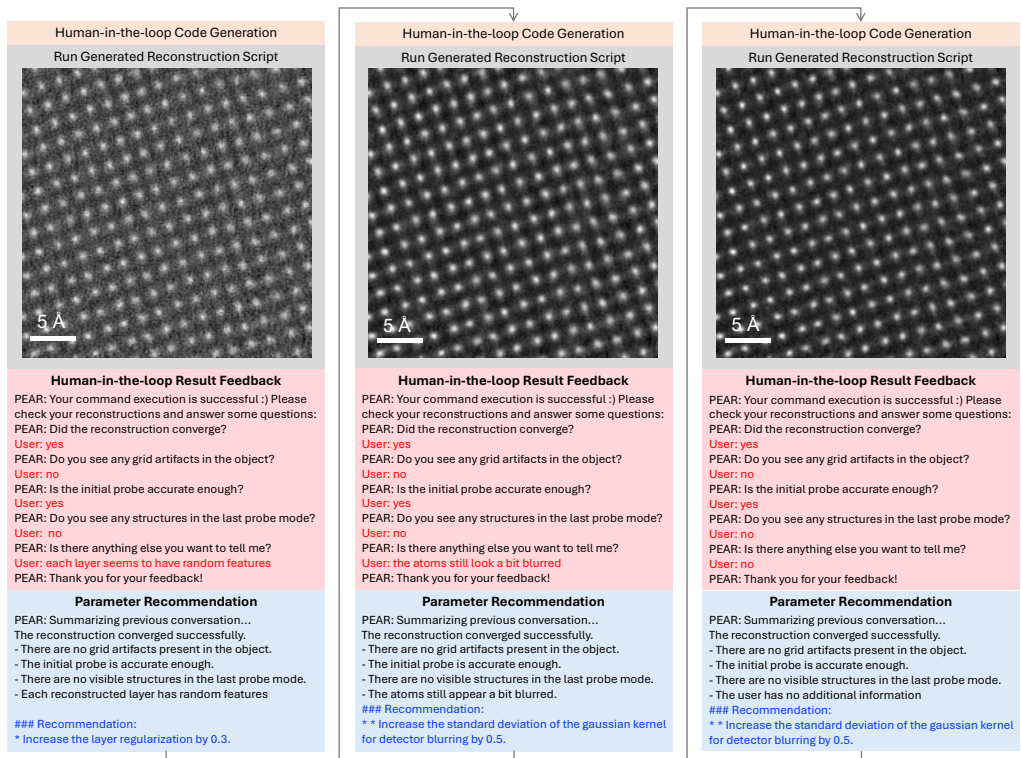


Figure 4: PEAR-assisted reconstructions of an electron ptychography dataset of SnSe. At automation Level 1, PEAR incorporates user feedback to recommend optimized parameters that enhance image quality. As a result, the reconstructed atomic structures become noticeably sharper and contain less noisy artifacts throughout the PEAR-guided workflow.

References

- Tomas Aidukas, Nicholas W Phillips, Ana Diaz, Emiliya Poghosyan, Elisabeth Müller, AFJ Levi, Gabriel Aeppli, Manuel Guizar-Sicairos, and Mirko Holler. High-performance 4-nm-resolution x-ray tomography using burst ptychography. *Nature*, 632(8023):81–88, 2024.
- Daniil A Boiko, Robert MacKnight, Ben Kline, and Gabe Gomes. Autonomous chemical research with large language models. *Nature*, 624(7992):570–578, 2023.
- Michael C Cao, Zhen Chen, Yi Jiang, and Yimo Han. Automatic parameter selection for electron ptychography via bayesian optimization. *Scientific Reports*, 12(1):12284, 2022.
- Zhen Chen, Michal Odstrcil, Yi Jiang, Yimo Han, Ming-Hui Chiu, Lain-Jong Li, and David A Muller. Dataset: Mixed-state electron ptychography enables sub-angstrom resolution imaging with picometer precision at low dose. *Nature Communications*, 2020.
- Matthew Cherukara, Alkaterini Vriza, Henry Chan, Tao Zhou, Varuni Katti Sastry, and Michael Prince. Calms: Context-aware language model for science. [Computer Software] <https://doi.org/10.11578/dc.20240410.1>, apr 2024. URL <https://doi.org/10.11578/dc.20240410.1>.
- Zehao Dong, Mengwu Huo, Jie Li, Jingyuan Li, Pengcheng Li, Hualei Sun, Lin Gu, Yi Lu, Meng Wang, Yayu Wang, et al. Visualization of oxygen vacancies and self-doped ligand holes in $\text{La}_3\text{Ni}_2\text{O}_7-\delta$. *Nature*, pages 1–6, 2024.
- Walter Hoppe. Beugung im inhomogenen primärstrahlwellenfeld. i. prinzip einer phasenmessung von elektronenbeugungsinterferenzen. *Acta Crystallographica Section A: Crystal Physics, Diffraction, Theoretical and General Crystallography*, 25(4):495–501, 1969.

- Yi Jiang, Zhen Chen, Yimo Han, Pratiti Deb, Hui Gao, Saien Xie, Prafull Purohit, Mark W Tate, Jiwoong Park, Sol M Gruner, et al. Electron ptychography of 2d materials to deep sub-ångström resolution. *Nature*, 559 (7714):343–349, 2018.
- Qiang Liu, Yu-Tong Liu, Chen Zhao, Qing-Song Weng, Junjing Deng, Inhui Hwang, Yi Jiang, Chengjun Sun, Tianyi Li, Wenqian Xu, et al. Conformal pedot coating enables ultra-high-voltage and high-temperature operation for single-crystal ni-rich cathodes. *ACS nano*, 16(9):14527–14538, 2022.
- Andres M. Bran, Sam Cox, Oliver Schilter, Carlo Baldassari, Andrew D White, and Philippe Schwaller. Augmenting large language models with chemistry tools. *Nature Machine Intelligence*, pages 1–11, 2024.
- John M Rodenburg and Helen ML Faulkner. A phase retrieval algorithm for shifting illumination. *Applied physics letters*, 85(20):4795–4797, 2004.
- Jiamei Song, Christopher S Allen, Si Gao, Chen Huang, Hidetaka Sawada, Xiaoqing Pan, Jamie Warner, Peng Wang, and Angus I Kirkland. Atomic resolution defocused electron ptychography at low dose with a fast, direct electron detector. *Scientific reports*, 9(1):3919, 2019.
- Dasol Yoon and David A Muller. Improving robustness of electron ptychography by bayesian optimization of tilt and thickness. *Microscopy and Microanalysis*, 30(Supplement_1), 2024.
- Dawei Zhang, Peggy Schoenherr, Pankaj Sharma, and Jan Seidel. Ferroelectric order in van der waals layered materials. *Nature Reviews Materials*, 8(1):25–40, 2023.

A Computational Experiment Setup

In the computational experiment described in Section 2.3, we used an advanced LLM to simulate human responses to the LLM agents’ questions, enabling us to efficiently conduct a large number of tests. Specifically, we employed the GPT-4 model via Azure OpenAI Service to simulate human responses. For the LLM agents, we utilized a range of models, including LLaMA 3.1 (8b, 70b, 405b) through Together AI, and ChatGPT models (gpt-4o, chatgpt-4o-latest, gpt-4o-mini) via OpenAI. The exact prompts used in these experiments are detailed in Figure 5.

The hypothetical experiment conditions and computational experiment logs are available in the supplementary materials.



Figure 5: Prompts used in the computational experiment.

B User Case Study Log

```
20240913-00:02:56 - Initializaing PEAR...
20240913-00:02:56 - ##### Settings #####
20240913-00:02:56 - Automation level: 1
20240913-00:02:56 - LLM: gpt-4o-mini
20240913-00:02:56 - Knowledge base: neurips_demo
20240913-00:02:56 - Data's base directory: /Modified/to/Hide/User/Info/
20240913-00:02:56 - Script directory: /Modified/to/Hide/User/Info/ptycho
20240913-00:02:56 - User name: User
20240913-00:02:56 - Computer name: lamda
20240913-00:02:56 - Matlab directory: /usr/local/bin/matlab
20240913-00:02:56 - External reconstruction script:
20240913-00:02:56 - #####
20240913-00:02:56 - PEAR: A driver file exists in the current working directory, it will be overwritten
.
20240913-00:02:56 - PEAR: Driver file created at driver.m
20240913-00:02:56 - PEAR: Hello User. Thank you for letting me assist with your ptychographic
reconstruction today.
20240913-00:02:56 - Agent: ParamsCollector
20240913-00:02:56 - PEAR: First I'd like to ask some questions about your data:
20240913-00:02:56 - PEAR: What is the scan number?
20240913-00:02:59 - User: 31
20240913-00:03:00 - PEAR: What is the beam energy (in keV)?
20240913-00:03:03 - User: 300
20240913-00:03:03 - PEAR: What is the radius of the bright field disk (in pixels)?
20240913-00:03:14 - User: 23.25
20240913-00:03:14 - PEAR: What is the convergence angle (in mrad)?
20240913-00:03:17 - User: 25
20240913-00:03:18 - PEAR: What is the detector size (in # of pixels)?
20240913-00:03:19 - User: 128
20240913-00:03:20 - PEAR: Do you want to load the initial object from an existing reconstruction? If so
, please provide the path to the reconstruction file.
20240913-00:03:22 - User: no
20240913-00:03:22 - PEAR: Do you want to generate the initial probe based on the ideal model or load it
from an existing reconstruction? If use ideal model, please provide the value of defocus (in
angstroms). If use existing probe, please provide the path to the reconstruction file.
20240913-00:03:44 - User: load an existing probe from /Modified/to/Hide/User/Info/Niter1000.mat
20240913-00:03:44 - PEAR: Is the initial probe accurate (similar size to the real probe)?
20240913-00:03:46 - User: yes
20240913-00:03:47 - PEAR: Do you want to generate scan positions assuming a perfect grid or load them
from an existing reconstruction? If you want to generate scan positions, please provide the scan
step size (in angstroms) and the number of scan points in both X and Y directions. If use existing
reconstructions, please provide the path to the reconstruction file.
20240913-00:04:05 - User: 128 x 128, step size is 0.2546 A
20240913-00:04:06 - PEAR: How many mixed-state probe modes do you want to use?
20240913-00:04:07 - User: 3
20240913-00:04:07 - PEAR: Which GPU do you want to use for reconstruction?
20240913-00:04:09 - User: 1
20240913-00:04:10 - PEAR: How many diffraction patterns does the scan have in total?
20240913-00:04:14 - User: 128 x 128
20240913-00:04:15 - PEAR: Did the sample drift during scan?
20240913-00:04:17 - User: no
20240913-00:04:17 - PEAR: What's the estimated sample thickness (in angstroms)?
20240913-00:04:20 - User: 100
20240913-00:04:20 - PEAR: Thank you for answering my questions.
20240913-00:04:20 - Agent: ParamsRecommender
20240913-00:04:20 - PEAR: Next I will initialize some reconstruction parameters based on your inputs.
20240913-00:04:20 - Agent: _summarize_conversation
20240913-00:04:20 - PEAR: Summarizing previous conversation...
20240913-00:04:23 - - The scan number is 31.
- The beam energy is 300 keV.
- The radius of the bright field disk is 23.25 pixels.
- The convergence angle is 25 mrad.
- The detector size is 128 pixels.
- The initial object will not be loaded from an existing reconstruction.
- An existing probe will be loaded from the path provided.
- The initial probe is accurate and similar in size to the real probe.
- The scan positions will be generated assuming a perfect grid with a step size of 0.2546 angstrom and
128 x 128 scan points.
- Three mixed-state probe modes will be used.
- GPU 1 will be used for reconstruction.
- The total number of diffraction patterns for the scan is 128 x 128.
- The sample did not drift during the scan.
- The estimated sample thickness is 100 angstrom.
20240913-00:04:25 - PEAR: Based on the provided information and the recommendations, here are the
suggested reconstruction parameters:

1. **Update Batch Size**: Set the update batch size to  $\lfloor \sqrt{128 \times 128} \rfloor = 128$ 
\).

2. **Gaussian Kernel Standard Deviation for Detector Blurring**: Since the beam energy is 300 keV, set
the standard deviation of the Gaussian kernel for detector blurring to **1**.

3. **Scan Position Correction**: Since the sample did not drift during the scan, set scan position
correction to **false**.

4. **Multisllice Ptychography**: Given that the estimated sample thickness is 100 angstrom, set
multisllice ptychography to **true** and the number of layers to **10** (1/10 of the sample
thickness).
```

```

5. **Number of Iterations**: Since the initial probe is accurate, set the number of iterations to
**50**.

These parameters should help optimize your ptychography reconstruction process.
20240913-00:04:25 - Agent: ParamsFormatter
20240913-00:04:28 - PEAR: Here are the current parameters:
20240913-00:04:28 - {
  "data_directory": "/Modified/to/Hide/User/Info/",
  "scan_number": 31,
  "beam_energy": 300,
  "radius_bright_field": 23.25,
  "convergence_angle": 25,
  "size_of_diffraction_patterns": 128,
  "use_external_object": false,
  "initial_object_path": "",
  "use_external_probe": true,
  "initial_probe_file": "/Modified/to/Hide/User/Info/Niter1000.mat",
  "defocus": 0,
  "use_external_positions": false,
  "initial_position_file": "",
  "grid_scan_positions": true,
  "scan_step_size_x": 0.2546,
  "scan_step_size_y": 0.2546,
  "number_scan_points_x": 128,
  "number_scan_points_y": 128,
  "number_of_iterations": 50,
  "update_batch_size": 128,
  "number_of_probe_modes": 3,
  "position_correction": false,
  "multislice_ptycho": true,
  "object_thickness": 100,
  "number_of_layers": 10,
  "layer_regularization_coefficient": 0,
  "diff_pattern_blur": 1,
  "gpu_id": 1
}
20240913-00:04:28 - ##### Reconstruction No.1: #####
20240913-00:04:28 - Agent: ParamsConfirmer
20240913-00:04:28 - PEAR: Let me know any changes you would like to make. If it looks good to you,
please say "LGTM".
20240913-00:04:51 - User: lgtm
20240913-00:04:51 - PEAR: Thank you for the feedback. I will create a reconstruction script using these
parameters.
20240913-00:04:51 - Agent: ScriptGenerator
20240913-00:04:51 - PEAR: Generating a reconstruction script...
20240913-00:04:57 - PEAR: The reconstruction script has been generated at /Modified/to/Hide/User/Info/
ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_1.m
20240913-00:04:57 - Agent: ScriptRunner
20240913-00:04:57 - PEAR: Running ptychographic reconstruction...
20240913-00:04:57 - PEAR: Executing script using /usr/local/bin/matlab on Modified-to-Hide-User-Info...
20240913-00:04:57 - PEAR: The MATLAB command used for reconstruction: /usr/local/bin/matlab -batch "
driver('/Modified/to/Hide/User/Info/ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_1.m')"
20240913-00:04:57 - PEAR: Executing command /usr/local/bin/matlab -batch "driver('/Modified/to/Hide/
User/Info/ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_1.m')"
20240913-00:12:23 - [init] : Preparing paths.
[Warning: Applying custom scan flip: 1 1 1 ]
[> In scans.positions.matlab_pos (line 14)
In scans.read_positions (line 20)
In core.initialize_ptycho (line 107)
In core.ptycho_recons/ptycho_call (line 172)
In core.ptycho_recons (line 289)
In ptycho_recon (line 576)
In ptycho_recon_pear (line 114)
In neurips_demo_gpt4o_mini_l1_ptycho_script_1 (line 44)
In run (line 99)
In driver (line 11)]
[init] : Preparing initial guess.
[init] : Preparing data using matlab APS data preparation.
[init] : Finished data preparation and initialization.
[ptycho] : Reconstructing S00031
[ptycho] : Calling engine GPU_MS
[GPU/CPU_MS-engine] : Intermediate results will be saved in the directory below every 50 iterations.
[GPU/CPU_MS-engine] : /Modified/to/Hide/User/Info/
[Warning: Applying custom data flip: 0 1 0 ]
[ptycho] : Elapsed time for engine GPU_MS: 425.7 s
Elapsed time is 438.322446 seconds.

20240913-00:12:23 - PEAR: Your command execution is successful :) Please take a look and give me some
feedback.
20240913-00:12:23 - Agent: QualityCollector
20240913-00:12:23 - PEAR: Please check your reconstructions and answer some questions:
20240913-00:12:24 - PEAR: Did the reconstruction converge?
20240913-00:12:44 - User: yes
20240913-00:12:44 - PEAR: Do you see any grid artifacts in the object?
20240913-00:12:47 - User: no
20240913-00:12:48 - PEAR: Is the initial probe accurate enough?
20240913-00:12:50 - User: yes
20240913-00:12:50 - PEAR: Do you see any structures in the last probe mode?
20240913-00:12:53 - User: yes
20240913-00:12:53 - PEAR: Is there anything else you want to tell me about the results?
20240913-00:13:01 - User: nothing at this moment
20240913-00:13:01 - PEAR: Thank you for your feedback!

```

```

20240913-00:13:01 - Agent: UpdatesRecommender
20240913-00:13:01 - PEAR: Next I will suggest some changes to the parameters:
20240913-00:13:01 - Agent: _summarize_conversation
20240913-00:13:01 - PEAR: Summarizing previous conversation...
20240913-00:13:02 - - The reconstruction converged successfully.
- There are no grid artifacts present in the object.
- The initial probe is accurate enough.
- Structures are visible in the last probe mode.
- There is nothing else to report about the results at this moment.
20240913-00:13:03 - Based on the CONVERSATION SUMMARY, the situation is as follows:

- The initial probe is accurate enough.
- The reconstruction has converged.
- The last probe mode has structures.

### Recommendation:
* Increase the number of probe modes by 3.
20240913-00:13:03 - Agent: ParamsUpdater
20240913-00:13:07 - PEAR: Here are the updated parameters that will be used to create the next
reconstruction script
20240913-00:13:07 - {
  "data_directory": "/Modified/to/Hide/User/Info/",
  "scan_number": 31,
  "beam_energy": 300,
  "radius_bright_field": 23.25,
  "convergence_angle": 25,
  "size_of_diffraction_patterns": 128,
  "use_external_object": false,
  "initial_object_path": "",
  "use_external_probe": true,
  "initial_probe_file": "/Modified/to/Hide/User/Info/Niter1000.mat",
  "defocus": 0,
  "use_external_positions": false,
  "initial_position_file": "",
  "grid_scan_positions": true,
  "scan_step_size_x": 0.2546,
  "scan_step_size_y": 0.2546,
  "number_scan_points_x": 128,
  "number_scan_points_y": 128,
  "number_of_iterations": 50,
  "update_batch_size": 128,
  "number_of_probe_modes": 6,
  "position_correction": false,
  "multisllice_ptycho": true,
  "object_thickness": 100,
  "number_of_layers": 10,
  "layer_regularization_coefficient": 0,
  "diff_pattern_blur": 1,
  "gpu_id": 1
}
20240913-00:13:07 - PEAR: Do you want to run another reconstruction? If so, please say "yes".
20240913-00:13:18 - User: yes
20240913-00:13:18 - ##### Reconstruction No.2: #####
20240913-00:13:18 - Agent: ParamsConfirmer
20240913-00:13:18 - PEAR: Let me know any changes you would like to make. If it looks good to you,
please say "LGTM".
20240913-00:13:21 - User: lgtm
20240913-00:13:21 - PEAR: Thank you for the feedback. I will create a reconstruction script using these
parameters.
20240913-00:13:21 - Agent: ScriptGenerator
20240913-00:13:21 - PEAR: Generating a reconstruction script...
20240913-00:13:26 - PEAR: The reconstruction script has been generated at /Modified/to/Hide/User/Info/
ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_2.m
20240913-00:13:26 - Agent: ScriptRunner
20240913-00:13:26 - PEAR: Running ptychographic reconstruction...
20240913-00:13:26 - PEAR: Executing script using /usr/local/bin/matlab on Modified-to-Hide-User-Info...
20240913-00:13:26 - PEAR: The MATLAB command used for reconstruction: /usr/local/bin/matlab -batch "
driver('/Modified/to/Hide/User/Info/ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_2.m')"
20240913-00:13:26 - PEAR: Executing command /usr/local/bin/matlab -batch "driver('/Modified/to/Hide/
User/Info/ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_2.m')"
20240913-00:25:12 - [init] : Preparing paths.
[Warning: Applying custom scan flip: 1 1 1 ]
[> In scans.positions.matlab_pos (line 14)
In scans.read_positions (line 20)
In core.initialize_ptycho (line 107)
In core.ptycho_recons/ptycho_call (line 172)
In core.ptycho_recons (line 289)
In ptycho_recon (line 576)
In ptycho_recon_pear (line 114)
In neurips_demo_gpt4o_mini_l1_ptycho_script_2 (line 44)
In run (line 99)
In driver (line 11)]
[init] : Preparing initial guess.
[init] : Preparing data using matlab APS data preparation.
[init] : Finished data preparation and initialization.
[ptycho] : Reconstructing S00031
[ptycho] : Calling engine GPU_MS
[GPU/CPU_MS-engine] : Intermediate results will be saved in the directory below every 50 iterations.
[GPU/CPU_MS-engine] : /Modified/to/Hide/User/Info/
[Warning: Applying custom data flip: 0 1 0 ]
[ptycho] : Elapsed time for engine GPU_MS: 684.4 s
Elapsed time is 697.096650 seconds.

```

```

20240913-00:25:12 - PEAR: Your command execution is successful :) Please take a look and give me some
feedback.
20240913-00:25:12 - Agent: QualityCollector
20240913-00:25:12 - PEAR: Please check your reconstructions and answer some questions:
20240913-00:25:13 - PEAR: Did the reconstruction converge?
20240913-00:25:27 - User: yes
20240913-00:25:27 - PEAR: Do you see any grid artifacts in the object?
20240913-00:25:30 - User: no
20240913-00:25:30 - PEAR: Is the initial probe accurate enough?
20240913-00:25:34 - User: yes
20240913-00:25:34 - PEAR: Do you see any structures in the last probe mode?
20240913-00:25:42 - User: not anymore
20240913-00:25:43 - PEAR: Is there anything else you want to tell me about the results?
20240913-00:26:19 - User: each layer in the reconstructed object seems to have random features that are
not real
20240913-00:26:20 - PEAR: Thank you for your feedback!
20240913-00:26:20 - Agent: UpdatesRecommender
20240913-00:26:20 - PEAR: Next I will suggest some changes to the parameters:
20240913-00:26:20 - Agent: _summarize_conversation
20240913-00:26:20 - PEAR: Summarizing previous conversation...
20240913-00:26:21 - - The reconstruction converged successfully.
- There are no grid artifacts present in the object.
- The initial probe is accurate enough.
- There are no visible structures in the last probe mode.
- Each layer in the reconstructed object seems to have random features that are not real.
20240913-00:26:21 - * Increase the layer regularization by 0.3.
20240913-00:26:21 - Agent: ParamsUpdater
20240913-00:26:25 - PEAR: Here are the updated parameters that will be used to create the next
reconstruction script
20240913-00:26:25 - {
  "data_directory": "/Modified/to/Hide/User/Info/",
  "scan_number": 31,
  "beam_energy": 300,
  "radius_bright_field": 23.25,
  "convergence_angle": 25,
  "size_of_diffraction_patterns": 128,
  "use_external_object": false,
  "initial_object_path": "",
  "use_external_probe": true,
  "initial_probe_file": "/Modified/to/Hide/User/Info/Niter1000.mat",
  "defocus": 0,
  "use_external_positions": false,
  "initial_position_file": "",
  "grid_scan_positions": true,
  "scan_step_size_x": 0.2546,
  "scan_step_size_y": 0.2546,
  "number_scan_points_x": 128,
  "number_scan_points_y": 128,
  "number_of_iterations": 100,
  "update_batch_size": 100,
  "number_of_probe_modes": 5,
  "position_correction": false,
  "multislice_ptycho": true,
  "object_thickness": 100,
  "number_of_layers": 10,
  "layer_regularization_coefficient": 0.3,
  "diff_pattern_blur": 1,
  "gpu_id": 1
}
20240913-00:26:25 - PEAR: Do you want to run another reconstruction? If so, please say "yes".
20240913-00:27:10 - User: yes
20240913-00:27:10 - ##### Reconstruction No.3: #####
20240913-00:27:10 - Agent: ParamsConfirmer
20240913-00:27:10 - PEAR: Let me know any changes you would like to make. If it looks good to you,
please say "LGTM".
20240913-00:27:23 - User: lgtm
20240913-00:27:23 - PEAR: Thank you for the feedback. I will create a reconstruction script using these
parameters.
20240913-00:27:23 - Agent: ScriptGenerator
20240913-00:27:23 - PEAR: Generating a reconstruction script...
20240913-00:27:29 - PEAR: The reconstruction script has been generated at /Modified/to/Hide/User/Info/
ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_3.m
20240913-00:27:29 - Agent: ScriptRunner
20240913-00:27:29 - PEAR: Running ptychographic reconstruction...
20240913-00:27:29 - PEAR: Executing script using /usr/local/bin/matlab on Modified-to-Hide-User-Info...
20240913-00:27:29 - PEAR: The MATLAB command used for reconstruction: /usr/local/bin/matlab -batch "
driver('/Modified/to/Hide/User/Info/ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_3.m')"
20240913-00:27:29 - PEAR: Executing command /usr/local/bin/matlab -batch "driver('/Modified/to/Hide/
User/Info/ptycho/neurips_demo_gpt4o_mini_l1_ptycho_script_3.m')"
20240913-00:49:23 - [init] : Preparing paths.
[Warning: Applying custom scan flip: 1 1 1 ]
[> In scans.positions.matlab_pos (line 14)
In scans.read_positions (line 20)
In core.initialize_ptycho (line 107)
In core.ptycho_recons/ptycho_call (line 172)
In core.ptycho_recons (line 289)
In ptycho_recon (line 576)
In ptycho_recon_pear (line 114)
In neurips_demo_gpt4o_mini_l1_ptycho_script_3 (line 44)
In run (line 99)
In driver (line 11)]

```

```

[init] : Preparing initial guess.
[init] : Preparing data using matlab APS data preparation.
[init] : Finished data preparation and initialization.
[ptycho] : Reconstructing S00031
[ptycho] : Calling engine GPU_MS
[GPU/CPU_MS-engine] : Intermediate results will be saved in the directory below every 50 iterations.
[GPU/CPU_MS-engine] : /Modified/to/Hide/User/Info/
[Warning: Applying custom data flip: 0 1 0 ]
[ptycho] : Elapsed time for engine GPU_MS: 1293.5 s
Elapsed time is 1306.163743 seconds.

20240913-00:49:23 - PEAR: Your command execution is successful :) Please take a look and give me some
feedback.
20240913-00:49:23 - Agent: QualityCollector
20240913-00:49:23 - PEAR: Please check your reconstructions and answer some questions:
20240913-00:49:23 - PEAR: Did the reconstruction converge?
20240913-00:56:04 - User: yes
20240913-00:56:05 - PEAR: Do you see any grid artifacts in the object?
20240913-00:56:15 - User: no
20240913-00:56:15 - PEAR: Is the initial probe accurate enough?
20240913-00:56:18 - User: yes
20240913-00:56:18 - PEAR: Do you see any structures in the last probe mode?
20240913-00:56:28 - User: no
20240913-00:56:28 - PEAR: Is there anything else you want to tell me about the results?
20240913-00:56:42 - User: the atoms still look a bit blurred
20240913-00:56:43 - PEAR: Thank you for your feedback!
20240913-00:56:43 - Agent: UpdatesRecommender
20240913-00:56:43 - PEAR: Next I will suggest some changes to the parameters:
20240913-00:56:43 - Agent: _summarize_conversation
20240913-00:56:43 - PEAR: Summarizing previous conversation...
20240913-00:56:44 - - The reconstruction converged successfully.
- There are no grid artifacts present in the object.
- The initial probe is accurate enough.
- There are no visible structures in the last probe mode.
- The atoms still appear a bit blurred.
20240913-00:56:44 - * Increase the standard deviation of the gaussian kernel for detector blurring by
0.5.
20240913-00:56:44 - Agent: ParamsUpdater
20240913-00:56:47 - PEAR: Here are the updated parameters that will be used to create the next
reconstruction script
20240913-00:56:47 - {
  "data_directory": "/Modified/to/Hide/User/Info/",
  "scan_number": 31,
  "beam_energy": 300,
  "radius_bright_field": 23.25,
  "convergence_angle": 25,
  "size_of_diffraction_patterns": 128,
  "use_external_object": false,
  "initial_object_path": "",
  "use_external_probe": true,
  "initial_probe_file": "/Modified/to/Hide/User/Info/Niter1000.mat",
  "defocus": 0,
  "use_external_positions": false,
  "initial_position_file": "",
  "grid_scan_positions": true,
  "scan_step_size_x": 0.2546,
  "scan_step_size_y": 0.2546,
  "number_scan_points_x": 128,
  "number_scan_points_y": 128,
  "number_of_iterations": 100,
  "update_batch_size": 100,
  "number_of_probe_modes": 5,
  "position_correction": false,
  "multisllice_ptycho": true,
  "object_thickness": 100,
  "number_of_layers": 10,
  "layer_regularization_coefficient": 0.3,
  "diff_pattern_blur": 1.5,
  "gpu_id": 1
}
20240913-00:56:47 - PEAR: Do you want to run another reconstruction? If so, please say "yes".
20240913-00:57:20 - User: yes
20240913-00:57:20 - ##### Reconstruction No.4: #####
20240913-00:57:20 - Agent: ParamsConfirmer
20240913-00:57:20 - PEAR: Let me know any changes you would like to make. If it looks good to you,
please say "LGTM".
20240913-00:57:34 - User: set the update batch size to 256
20240913-00:57:38 - PEAR: Here are the updated parameters:
20240913-00:57:38 - {
  "data_directory": "/Modified/to/Hide/User/Info/",
  "scan_number": 31,
  "beam_energy": 300,
  "radius_bright_field": 23.25,
  "convergence_angle": 25,
  "size_of_diffraction_patterns": 128,
  "use_external_object": false,
  "initial_object_path": "",
  "use_external_probe": true,
  "initial_probe_file": "/Modified/to/Hide/User/Info/Niter1000.mat",
  "defocus": 0,
  "use_external_positions": false,
  "initial_position_file": "",

```



```

    "grid_scan_positions": true,
    "scan_step_size_x": 0.2546,
    "scan_step_size_y": 0.2546,
    "number_scan_points_x": 128,
    "number_scan_points_y": 128,
    "number_of_iterations": 100,
    "update_batch_size": 256,
    "number_of_probe_modes": 5,
    "position_correction": false,
    "multislice_ptycho": true,
    "object_thickness": 100,
    "number_of_layers": 10,
    "layer_regularization_coefficient": 0.3,
    "diff_pattern_blur": 1.5,
    "gpu_id": 1
}
20240913-00:57:38 - PEAR: Let me know any changes you would like to make. If it looks good to you,
    please say "LGTM".
20240913-00:57:50 - User: lgtm
20240913-00:57:50 - PEAR: Thank you for the feedback. I will create a reconstruction script using these
    parameters.
20240913-00:57:50 - Agent: ScriptGenerator
20240913-00:57:50 - PEAR: Generating a reconstruction script...
20240913-00:57:56 - PEAR: The reconstruction script has been generated at /Modified/to/Hide/User/Info/
    ptycho/neurips_demo_gpt4o_mini_ll_ptycho_script_4.m
20240913-00:57:56 - Agent: ScriptRunner
20240913-00:57:56 - PEAR: Running ptychographic reconstruction...
20240913-00:57:56 - PEAR: Executing script using /usr/local/bin/matlab on Modified-to-Hide-User-Info...
20240913-00:57:56 - PEAR: The MATLAB command used for reconstruction: /usr/local/bin/matlab -batch "
    driver('/Modified/to/Hide/User/Info/ptycho/neurips_demo_gpt4o_mini_ll_ptycho_script_4.m')"
20240913-00:57:56 - PEAR: Executing command /usr/local/bin/matlab -batch "driver('/Modified/to/Hide/
    User/Info/ptycho/neurips_demo_gpt4o_mini_ll_ptycho_script_4.m')"
20240913-01:17:09 - [init] : Preparing paths.
[Warning: Applying custom scan flip: 1 1 1 ]
[> In scans.positions.matlab_pos (line 14)
In scans.read_positions (line 20)
In core.initialize_ptycho (line 107)
In core.ptycho_recons/ptycho_call (line 172)
In core.ptycho_recons (line 289)
In ptycho_recon (line 576)
In ptycho_recon_pear (line 114)
In neurips_demo_gpt4o_mini_ll_ptycho_script_4 (line 44)
In run (line 99)
In driver (line 11)]
[init] : Preparing initial guess.
[init] : Preparing data using matlab APS data preparation.
[init] : Finished data preparation and initialization.
[ptycho] : Reconstructing S00031
[ptycho] : Calling engine GPU_MS
[GPU/CPU_MS-engine] : Intermediate results will be saved in the directory below every 50 iterations.
[GPU/CPU_MS-engine] : /Modified/to/Hide/User/Info/
[Warning: Applying custom data flip: 0 1 0 ]
[ptycho] : Elapsed time for engine GPU_MS: 1132.4 s
Elapsed time is 1145.022578 seconds.

20240913-01:17:09 - PEAR: Your command execution is successful :) Please take a look and give me some
    feedback.
20240913-01:17:09 - Agent: QualityCollector
20240913-01:17:09 - PEAR: Please check your reconstructions and answer some questions:
20240913-01:17:09 - PEAR: Did the reconstruction converge?
20240913-01:19:28 - User: yes
20240913-01:19:28 - PEAR: Do you see any grid artifacts in the object?
20240913-01:19:30 - User: no
20240913-01:19:31 - PEAR: Is the initial probe accurate enough?
20240913-01:19:36 - User: yes
20240913-01:19:36 - PEAR: Do you see any structures in the last probe mode?
20240913-01:19:39 - User: no
20240913-01:19:40 - PEAR: Is there anything else you want to tell me about the results?
20240913-01:19:42 - User: no
20240913-01:19:43 - PEAR: Thank you for your feedback!
20240913-01:19:43 - Agent: UpdatesRecommender
20240913-01:19:43 - PEAR: Next I will suggest some changes to the parameters:
20240913-01:19:43 - Agent: _summarize_conversation
20240913-01:19:43 - PEAR: Summarizing previous conversation...
20240913-01:19:44 - - The reconstruction converged successfully.
- There are no grid artifacts present in the object.
- The initial probe is accurate enough.
- There are no structures visible in the last probe mode.
- The user has no additional information to provide about the results.
- All diagnosis questions have been asked.
20240913-01:19:45 - * Increase the standard deviation of the gaussian kernel for detector blurring by
    0.5.
20240913-01:19:45 - Agent: ParamsUpdater
20240913-01:19:48 - PEAR: Here are the updated parameters that will be used to create the next
    reconstruction script
20240913-01:19:48 - {
    "data_directory": "/Modified/to/Hide/User/Info/",
    "scan_number": 31,
    "beam_energy": 300,
    "radius_bright_field": 23.25,
    "convergence_angle": 25,
    "size_of_diffraction_patterns": 128,

```

```
"use_external_object": false,
"initial_object_path": "",
"use_external_probe": true,
"initial_probe_file": "/Modified/to/Hide/User/Info//Niter1000.mat",
"defocus": 0,
"use_external_positions": false,
"initial_position_file": "",
"grid_scan_positions": true,
"scan_step_size_x": 0.2546,
"scan_step_size_y": 0.2546,
"number_scan_points_x": 128,
"number_scan_points_y": 128,
"number_of_iterations": 100,
"update_batch_size": 256,
"number_of_probe_modes": 5,
"position_correction": false,
"multislice_ptycho": true,
"object_thickness": 100,
"number_of_layers": 10,
"layer_regularization_coefficient": 0.3,
"diff_pattern_blur": 2.0,
"gpu_id": 1
}
20240913-01:19:48 - PEAR: Do you want to run another reconstruction? If so, please say "yes".
20240913-01:20:20 - User: no
20240913-01:20:20 - PEAR: Reconstructions completed. Have a nice day!
```