

# Revolutionizing Cross-sectional Imaging — Supplementary Materials

Y. F. Fan<sup>1,3\*</sup>, Y. B. Fan<sup>2,3\*</sup>, L. P. Luo<sup>4\*</sup>, W. T. Lin<sup>1</sup>, Z. Y. Li<sup>5</sup>, X.

Zhong<sup>4</sup>, C. Z. Shi<sup>4</sup>, T. Newman<sup>5</sup>, Y. Zhou<sup>1</sup>, C. S. Lv<sup>1</sup> and Y. Z. Fan<sup>5</sup>

<sup>1</sup>Center for Scientific Research, Guangzhou Institute of Physical Education, Guangzhou 510500, P.R. China

<sup>2</sup>Bioengineering Department, Beijing University of Aeronautics and Astronautics, Beijing 100191, P.R. China

<sup>3</sup>Key Laboratory of Optimal Design and Evaluation of Medical Equipment and Rehabilitation Aids, BUAA Research Institute, Guangzhou 510530, P.R. China

<sup>4</sup>Medical Imaging Center, the First Affiliated Hospital of Jinan University, Guangzhou 510632, P.R. China

<sup>5</sup>College of Foreign Studies, Jinan University, Guangzhou 510632, P.R. China

(Dated: December 3, 2024)

PACS numbers: 87.85.Pq, 42.30.Wb.

## I. CONVERTING CROSS-SECTIONAL IMAGE INTO TEXT FORMAT

3D Medical Image Processing and Analyzing System (3DMed) [1] is a free software. Mimics may convert better, but it is not free.

## II. CALCULATING THE CENTER OF MASS (COM)

In Worksheets "Sesamoid (left foot)" and "Sesamoid (right foot)", COM is calculated by the following equation:

$$\left( x_c = \frac{\sum_{i=1}^n (x_i g_i)}{\sum_{i=1}^n g_i}, y_c = \frac{\sum_{i=1}^n (y_i g_i)}{\sum_{i=1}^n g_i}, z_c = \frac{\sum_{i=1}^n (z_i g_i)}{\sum_{i=1}^n g_i} \right),$$

where  $(x_c, y_c, z_c)$  refers to COM position,  $(x_i, y_i, z_i)$  the position of pixel unit,  $g_i$  the gray value of pixel unit and  $n$  the count of pixel units. Please refer to Literature [2, 3] for the relationship between the gray value and its density.

## III. POSITIONING

In appendix *Software S1.xlsx*, we copied range A:D of Worksheet "Sesamoid (left foot)" or Worksheet "Sesamoid (right foot)" to range A:D of Worksheet "Data" to calculate the position in Worksheet "Position".

### A. Method

In Worksheet "Position", cells (A1,B1,C1) refer to the COM position (unit: mm); cell D1 the mean value of gray values, cell D2 the count of pixel unit; cells A4, B4 the horizontal and vertical distance between pixel (unit: mm), cell C4 slices distance (unit: mm) and cell D3 the count +4 of pixel unit. Range A5:D(D3) includes data from Worksheet "Data", where range A5:A(D3): Coordinate  $x$  (mm); range B5:B(D3): Coordinate  $y$  (mm); range C5:C(D3): Coordinate  $z$  (mm); range D5:D(D3): Gray values. Range F5:F(D3): Coordinate  $y$  (mm); range G5:G(D3): Coordinate  $z$  (mm). They are the new coordinates of pixel units after rotating about axis  $x$ . Range H5:H(D3): Coordinate  $x$  (mm); range J5:J(D3): Coordinate  $z$  (mm). They are the new coordinates of pixel units after rotating about axis  $y$ . Range K5:K(D3): Coordinate  $x$  (mm); range L5:L(D3): Coordinate  $y$  (mm). They are the new coordinates of pixel units after rotating about axis  $z$ . The rotation about axis  $x$  is performed by the following equation:

$$\begin{pmatrix} x_i^1 \\ y_i^1 \\ z_i^1 \end{pmatrix} = \begin{pmatrix} 1 & 0 & 0 \\ 0 & \cos(E5) & -\sin(E5) \\ 0 & \sin(E5) & \cos(E5) \end{pmatrix} \begin{pmatrix} x_i^{data} \\ y_i^{data} - y_c \\ z_i^{data} - z_c \end{pmatrix} + \begin{pmatrix} 0 \\ y_c \\ z_c \end{pmatrix}, \quad (1)$$

\* Corresponding authors

Email: tfyf@gipe.edu.cn, yubofan@buaa.edu.cn and tluolp@jnu.edu.cn

where E5 stands for the angular displacement,  $(x_i^{data}, y_i^{data}, z_i^{data})$  the pixel unit coordinates in Worksheet "Data",  $(x_c, y_c, z_c)$  the position of COM, and  $(x_i^1, y_i^1, z_i^1)$  pixel unit coordinates after rotation about axis  $x$ .

The rotation about axis  $y$  is performed by the following equation:

$$\begin{pmatrix} x_i^2 \\ y_i^2 \\ z_i^2 \end{pmatrix} = \begin{pmatrix} \cos(F5) & 0 & \sin(F5) \\ 0 & 1 & 0 \\ -\sin(F5) & 0 & \cos(F5) \end{pmatrix} \begin{pmatrix} x_i^1 - x_c \\ y_i^1 \\ z_i^1 - z_c \end{pmatrix} + \begin{pmatrix} x_c \\ 0 \\ z_c \end{pmatrix}, \quad (2)$$

where F5 stands for the angular displacement.  $(x_c, y_c, z_c)$  and  $(x_i^1, y_i^1, z_i^1)$  have the same definitions as those in Eq. (1), and  $(x_i^2, y_i^2, z_i^2)$  pixel unit coordinates after rotation about axis  $y$ .

The rotation about axis  $z$  is performed by the following equation:

$$\begin{pmatrix} x_i^3 \\ y_i^3 \\ z_i^3 \end{pmatrix} = \begin{pmatrix} \cos(G5) & -\sin(G5) & 0 \\ \sin(G5) & \cos(G5) & 0 \\ 0 & 0 & 1 \end{pmatrix} \begin{pmatrix} x_i^2 - x_c \\ y_i^2 - y_c \\ z_i^2 \end{pmatrix} + \begin{pmatrix} x_c \\ y_c \\ 0 \end{pmatrix}, \quad (3)$$

where G5 stands for the angular displacement.  $(x_c, y_c, z_c)$ ,  $(x_i^2, y_i^2, z_i^2)$  have the same definitions as those in Eq. (2), and  $(x_i^3, y_i^3, z_i^3)$  pixel unit coordinates after rotation about axis  $z$ .

In Eq. (2) - (3),

$$\begin{cases} E5 = K5 = \frac{1}{2} \arctan\left(\frac{2G1}{(E1-F1)}\right) \\ F5 = L5 = \frac{1}{2} \arctan\left(\frac{2G2}{(E2-F2)}\right) \\ G5 = M5 = \frac{1}{2} \arctan\left(\frac{2G3}{(E3-F3)}\right) \end{cases}. \quad (4)$$

In Eq. (4),

$$\begin{cases} E1 = \sum_{i=1}^n (y_i^{data} - y_c)^2 \Delta V D_i \\ E2 = \sum_{i=1}^n (x_i^1 - x_c)^2 \Delta V D_i \\ E3 = \sum_{i=1}^n (x_i^2 - x_c)^2 \Delta V D_i \end{cases},$$

$$\begin{cases} F1 = \sum_{i=1}^n (z_i^{data} - z_c)^2 \Delta V D_i \\ F2 = \sum_{i=1}^n (z_i^1 - z_c)^2 \Delta V D_i \\ F3 = \sum_{i=1}^n (y_i^2 - y_c)^2 \Delta V D_i \end{cases},$$

$$\begin{cases} G1 = \sum_{i=1}^n (y_i^{data} - y_c) (z_i^{data} - z_c) \Delta V D_i \\ G2 = \sum_{i=1}^n (x_i^1 - x_c) (z_i^1 - z_c) \Delta V D_i \\ G3 = \sum_{i=1}^n (x_i^2 - x_c) (y_i^2 - y_c) \Delta V D_i \end{cases},$$

where  $\Delta V = A4B4C4$  (volume element).

Cell N1 stands for the moment of inertia (MOI) before rotation, cell N2 MOI after rotation about axis  $x$ , cell N3 MOI after rotation about axis  $y$ , cell N4 MOI after rotation about axis  $z$ . The MOI of sesamoid is calculated by the following equation:

$$\begin{cases} N1 = H1 + I1 + J1 \\ N2 = H2 + I2 + J2 \\ N3 = H3 + I3 + J3 \\ N4 = H4 + I4 + J4 \end{cases}, \quad (5)$$

where

$$\begin{cases} H1 = \sum_{i=1}^n \left( (y_i^{data} - y_c)^2 + (z_i^{data} - z_c)^2 \right) \Delta V D_i \\ H2 = \sum_{i=1}^n \left( (y_i^1 - y_c)^2 + (z_i^1 - z_c)^2 \right) \Delta V D_i \\ H3 = \sum_{i=1}^n \left( (y_i^2 - y_c)^2 + (z_i^2 - z_c)^2 \right) \Delta V D_i \\ H4 = \sum_{i=1}^n \left( (y_i^3 - y_c)^2 + (z_i^3 - z_c)^2 \right) \Delta V D_i \end{cases},$$

TABLE I: First positioning of Sesamoid (left foot)

	E	F	G	H	I	J	N
1	2.05	1.40	0.41	3.45	2.20	2.85	8.50
2	0.80	1.20	-0.17	3.45	2.00	3.05	8.50
3	0.74	2.25	0.07	3.51	2.00	2.99	8.50
4	-26.01	20.23	2.60	3.51	2.00	2.99	8.50

For range E1:G3, unit is  $g \cdot cm^2$ ; for cell E4, F4, G4, unit is degree; for range H1:N4, unit is  $g \cdot cm^2$ .

$$\left\{ \begin{array}{l} I1 = \sum_{i=1}^n \left( (x_i^{data} - x_c)^2 + (z_i^{data} - z_c)^2 \right) \Delta V D_i \\ I2 = \sum_{i=1}^n \left( (x_i^1 - x_c)^2 + (z_i^1 - z_c)^2 \right) \Delta V D_i \\ I3 = \sum_{i=1}^n \left( (x_i^2 - x_c)^2 + (z_i^2 - z_c)^2 \right) \Delta V D_i \\ I4 = \sum_{i=1}^n \left( (x_i^3 - x_c)^2 + (z_i^3 - z_c)^2 \right) \Delta V D_i \end{array} \right. ,$$

$$\left\{ \begin{array}{l} J1 = \sum_{i=1}^n \left( (x_i^{data} - x_c)^2 + (y_i^{data} - y_c)^2 \right) \Delta V D_i \\ J2 = \sum_{i=1}^n \left( (x_i^1 - x_c)^2 + (y_i^1 - y_c)^2 \right) \Delta V D_i \\ J3 = \sum_{i=1}^n \left( (x_i^2 - x_c)^2 + (y_i^2 - y_c)^2 \right) \Delta V D_i \\ J4 = \sum_{i=1}^n \left( (x_i^3 - x_c)^2 + (y_i^3 - y_c)^2 \right) \Delta V D_i \end{array} \right. ,$$

$(x_c, y_c, z_c)$ ,  $(x_i^1, y_i^1, z_i^1)$ ,  $(x_i^2, y_i^2, z_i^2)$  and  $(x_i^3, y_i^3, z_i^3)$  have the same definitions as those in Eq. (1) - (4).

## B. Operation and results

### 1. Positioning of sesamoid (left foot)

Select Worksheet "Sesamoid (left foot)" range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Data" range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Paste**.

Revise some parameters in Worksheet "Position". Select range 1:4, on the **Home tab**, in the **Editing group**, and then click **Find Select**, click **Replace**, in Find what: a certain number, and Replace with: 4524, click **Replace All**. The results are shown in Table (I):

Select Worksheet "Position" range K:M, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Temporary data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste Special**, click **Values** (Paste), click **None** (Operation), click **OK**. Select Worksheet "Temporary data" range 1:4, on the **Home tab**, in the **Cells group**, and then click **Delete**, click **Delete** Sheet Rows. Select Worksheet "Temporary data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste**. The results are shown in Table (II).

Repeat the operation of Table (II) and then get Table (III).

### 2. Positioning of sesamoid (right foot)

Repeat the above operations of Tables (I), (II) and (III) on sesamoid (right foot). We get Tables (IV), (V), and (VI) respectively.

TABLE II: Second positioning of sesamoid (left foot)

	E	F	G	H	I	J	N
1	2.25	1.26	-0.03	3.51	2.00	2.99	8.50
2	0.74	1.26	0.00	3.51	2.00	2.99	8.50
3	0.74	2.25	0.00	3.51	2.00	2.99	8.50
4	1.46	-0.12	-0.00	3.51	2.00	2.99	8.50

The unit of each item is the same as that in Table (I).

TABLE III: Third positioning of sesamoid (left foot)

	E	F	G	H	I	J	N
1	2.25	1.26	0.00	3.51	2.00	2.99	8.50
2	0.74	1.26	0.00	3.51	2.00	2.99	8.50
3	0.74	2.25	0.00	3.51	2.00	2.99	8.50
4	0.00	0.00	0.00	3.51	2.00	2.99	8.50

The unit of each item is the same as that in Table (I).

Revise some parameters in Worksheet "Position". Select range 1:4, on the **Home tab**, in the **Editing group**, and then click Find Select, click **Replace**, in Find what: 4524, and Replace with: 5044, click **Replace All**. The results are shown in Table (IV):

Fig. 1 shows the result after the sesamoid is positioned.

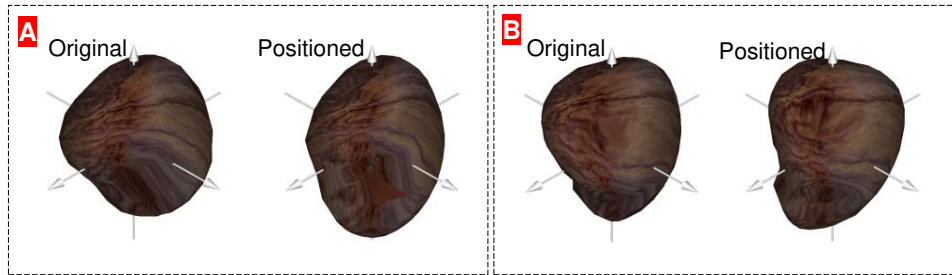


FIG. 1: Positioning of sesamoid. (A) Body coordinates of sesamoid (left foot) before and after being positioned; (B) Body coordinates of sesamoid (right foot) before and after being positioned.

TABLE IV: First positioning of sesamoid (right foot)

	E	F	G	H	I	J	N
1	2.58	1.76	0.46	4.34	2.70	3.52	10.56
2	0.94	1.55	0.16	4.34	2.49	3.73	10.56
3	0.90	2.79	0.05	4.38	2.49	3.69	10.56
4	-24.14	-13.71	1.66	4.38	2.49	3.69	10.56

The unit of each item is the same as that in Table (I).

TABLE V: Second positioning of sesamoid (right foot)

	E	F	G	H	I	J	N
1	2.79	1.59	0.01	4.38	2.49	3.69	10.56
2	0.90	1.59	0.00	4.38	2.49	3.69	10.56
3	0.90	2.79	0.00	4.38	2.49	3.69	10.56
4	-0.64	0.03	-0.00	4.38	2.49	3.69	10.56

The unit of each item is the same as that in Table (I).

TABLE VI: Third positioning of sesamoid (right foot)

	E	F	G	H	I	J	N
1	2.79	1.59	0.00	4.38	2.49	3.69	10.56
2	0.90	1.59	0.00	4.38	2.49	3.69	10.56
3	0.90	2.79	0.00	4.38	2.49	3.69	10.56
4	-0.00	-0.00	0.00	4.38	2.49	3.69	10.56

The unit of each item is the same as that in Table (I).

### 3. Positioning of sesamoid (right foot) after the coordinate transformation

To further demonstrate the function of positioning, let the positioned sesamoid (right foot) rotate about axis x, y, and z at 3, 5, and 7 degrees respectively. Observe the positioning change after the coordinate transformation.

In Worksheet "Position", let cell E4=3(degree), cell F4=5(degree), cell G4=7(degree).

Select Worksheet "Position" range K:M, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Temporary data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste Special**, click **Values** (Paste), click **None** (Operation), click **OK**.

Select Worksheet "Temporary data" range 1:4, on the **Home tab**, in the Cells group, and then click **Delete**, click **Delete** Sheet Rows. Select Worksheet "Temporary data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste**.

In Worksheet "Position", let cell E4=A3, cell F4=B3, cell G4=C3. Repeat the above operations of Tables (I), (II) and (III) on sesamoid (right foot). We get Tables (VII), (VIII), and (IX) respectively.

In Worksheet "Position" of appendix Sesamoid.xlsb, cells N1-N4 must be a constant because the MOI is irrespective

TABLE VII: First positioning of sesamoid (right foot)

	E	F	G	H	I	J	N
1	2.76	1.58	0.07	4.35	2.52	3.69	10.56
2	0.93	1.58	0.06	4.35	2.51	3.70	10.56
3	0.93	2.77	-0.22	4.35	2.51	3.69	10.56
4	-3.37	-5.65	-6.66	4.38	2.49	3.69	10.56

The unit of each item is the same as that in Table (I).

TABLE VIII: Second positioning of sesamoid (right foot)

	E	F	G	H	I	J	N
1	2.79	1.59	-0.02	4.38	2.49	3.69	10.56
2	0.90	1.59	0.00	4.38	2.49	3.69	10.56
3	0.90	2.79	0.00	4.38	2.49	3.69	10.56
4	1.02	0.21	0.00	4.38	2.49	3.69	10.56

The unit of each item is the same as that in Table (I).

TABLE IX: Third positioning of sesamoid (right foot)

	E	F	G	H	I	J	N
1	2.79	1.59	0.00	4.38	2.49	3.69	10.56
2	0.90	1.59	0.00	4.38	2.49	3.69	10.56
3	0.90	2.79	0.00	4.38	2.49	3.69	10.56
4	0.00	-0.00	-0.00	4.38	2.49	3.69	10.56

The unit of each item is the same as that in Table (I).

to the transformation of coordinates when MOI goes/is on tissue or organ. In addition, the positioning indices of cells H1=H2=H3, I1=I2=I3 and J1=J2=J3 will result in ells H1=H2=H3, I1=I2=I3 and J1=J2=J3, which serves as a condition whether it is positioned.

Please refer to Literature [4] for rotation and positioning.

#### IV. RECONSTRUCTION

In Worksheet "Reconstruction", range A-D includes the coordinates and grey values of the pixel units. (I1,J1,K1) refers to the COMs of sesamoid while (I2,J2,K2) the COMs of reconstructed sesamoid.

The reconstructed CSI is isotropic:

$$\begin{cases} x'_i = \text{Round} \left( \left( \frac{(x_i - x_c)}{R}, 0 \right) R \right) \\ y'_i = \text{Round} \left( \left( \frac{(y_i - y_c)}{R}, 0 \right) R \right) \\ z'_i = \text{Round} \left( \left( \frac{(z_i - z_c)}{R}, 0 \right) R \right) \end{cases}, \quad (6)$$

where  $R$  stands for CSI resolution (0.5mm).

The following is a cross section along one axis.

Along anteroposterior axis,

$$\begin{cases} x'_i = \begin{cases} x_i & p + \Delta > x_i > p - \Delta \\ \text{space} & \text{other} \end{cases} \\ y'_i = \begin{cases} y_i & \text{other} \\ \text{space} & x'_i = \text{space} \end{cases} \\ z'_i = \begin{cases} z_i & \text{other} \\ \text{space} & x'_i = \text{space} \end{cases} \end{cases}. \quad (7)$$

Along left-right axis,

$$\left\{ \begin{array}{l} x'_i = \begin{cases} x_i & \text{other} \\ space & y'_i = space \end{cases} \\ y'_i = \begin{cases} y_i & p + \Delta > y_i > p - \Delta \\ space & \text{other} \end{cases} \\ z'_i = \begin{cases} z_i & \text{other} \\ space & y'_i = space \end{cases} \end{array} \right. . \quad (8)$$

Along vertical axis,

$$\left\{ \begin{array}{l} x'_i = \begin{cases} x_i & \text{other} \\ space & z'_i = space \end{cases} \\ y'_i = \begin{cases} y_i & \text{other} \\ space & z'_i = space \end{cases} \\ z'_i = \begin{cases} z_i & p + \Delta > z_i > p - \Delta \\ space & \text{other} \end{cases} \end{array} \right. . \quad (9)$$

In Worksheets "Anteroposterior axis", "Left-right axis", and "Vertical axis", cell E1 stands for  $p$  in Eq. (7) - (9),  $p$  a point on axis,  $p \in (0, \pm R, \pm 2R \dots)$ .  $R$  is defined in (6), and  $0 < \Delta < \frac{1}{2}R$ . This is to reconstruct a section.

In Worksheet "Reconstruction", the values of cells I2, J2, and K2 reflect the effect of reconstruction. To minimize the value is a basic principle in reconstruction. In Worksheets "Anteroposterior axis", "Left-right axis", and "Vertical axis", I1 is essential because it determines whether the data retrieved are from the same section. In Worksheets "Coronal Map", "Sagittal map", and "Horizontal map", the recovering of the lost pixel unit needs further research.

In Worksheets "Coronal Section", "Sagittal Section" and "Horizontal Section", more research is needed to investigate the reconstruction; for example, the reconstruction of between two adjacent sections is not discussed in this appendix.

## V. RECONSTRUCTING THE INTERESTED CROSS SECTIONS

In this paper, region of interest (ROI) refers to the section along certain axis that a researcher or scientist wants to examine. It is arbitrary.

### A. COM on three sections

In Worksheet "Anteroposterior axis", cell E1=0. Then, Select range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Temporary data" range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste Special**, click **Values** (Paste), click **None** (Operation), click **OK**.

Select Worksheet "Temporary data" range B:B, on the **Home tab**, in the **Editing group**, and click Sort Filter, click Sort Smallest to Largest, click Expand the selection, click Sort. Select Worksheet "Temporary data" range C:C, on the **Home tab**, in the **Editing group**, and click Sort Filter, click Sort Smallest to Largest, click Expand the selection, click Sort.

Select Worksheet "Temporary data" range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Coronal Section" range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste**.

Select Worksheet "Coronal Section" range N:P, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Temporary data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste Special**, click **Values** (Paste), click **None** (Operation), click **OK**.

Select Worksheet "Temporary data" range A:A, on the **Home tab**, in the **Editing group**, and click Sort Filter, click Sort Smallest to Largest, click Expand the selection, click Sort. Select Worksheet "Temporary data" range B:B, on the **Home tab**, in the **Editing group**, and click Sort Filter, click Sort Smallest to Largest, click Expand the selection, click Sort.

Select Worksheet "Temporary data" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Three planes" range A:C, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste**.

With the same operation, we can reconstruct sections on the Sagittal Plane and Horizontal Plane. Worksheet "Three planes" respectively: range A:C are Coronal Section, range D:F Sagittal Section, and range G:I Horizontal Section.

This will successfully retrieve Coronal Section. In Worksheet "Coronal Section", range E:M can eliminate pixel units that have same coordinates.

For Worksheet "Coronal Map", retrieve Coronal Section from Worksheet "Three planes". In Worksheet "Coronal Map", Worksheet "Sagittal map" and Worksheet "Horizontal map", range F1:EY50, change y coordinates of the section pixel units to columns in Worksheet, range F51:EY100. Complete the section pixel units coordinates values and change them into integers.

Range F101:EY250. Retrieve x and y coordinate values of section pixel units and transform the row and column number. The cell value will be the grey value of section pixel units. Range H201:BE250, map the section; range H251:BE300, recover the lost points of section pixel units. Thus map the Coronal Section.

With the same method/procedure, map the other two sections from Worksheet "Sagittal map" and Worksheet "Horizontal map". The results are shown in Worksheet "Maps". See Fig. 2.

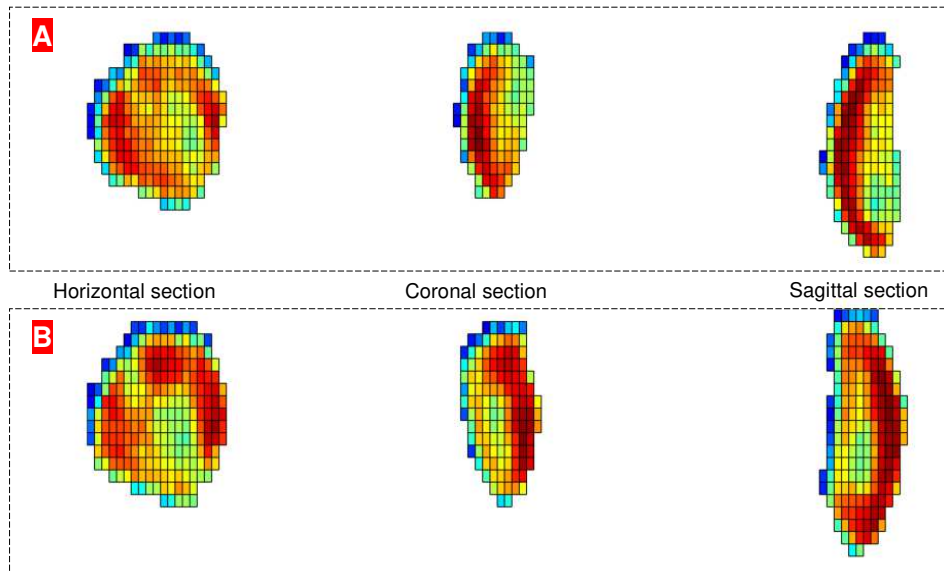


FIG. 2: Sesamoid section after being positioned. (A) Three sections of sesamoid (left foot). (B) Three sections of sesamoid (right foot).

### B. Any section along any axis

After Sesamoid (left foot) is positioned, select Worksheet "Result" range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Copy**. Select Worksheet "Data" range A:D, on the **Home tab**, in the **Clipboard group**, and then click **Paste**, click **Paste Special**, click **Values** (Paste), click **None** (Operation), click **OK**.

Revise some parameters in Worksheet "Position". Select range 1:4, on the **Home tab**, in the **Editing group**, and then click Find Select, click **Replace**, in Find what: a certain number, and Replace with: 4524, click **Replace All**.

In Worksheet "Position", let cell E4=12 (degree); in Worksheet "Anteroposterior axis", let cell E1=2 (mm); in Worksheet "Left-right axis" and Worksheet "Vertical axis", let cell E1=0 (mm). In this way, Anteroposterior axis has rotates 12 degrees. In the new coordinate system, retrieve a section from Coronal Section vertical to anteroposterior axis at a distance of 2mm from COM.

Set the same parameters, transform the coordinate of Sesamoid (right foot). See Fig. 3 for results.

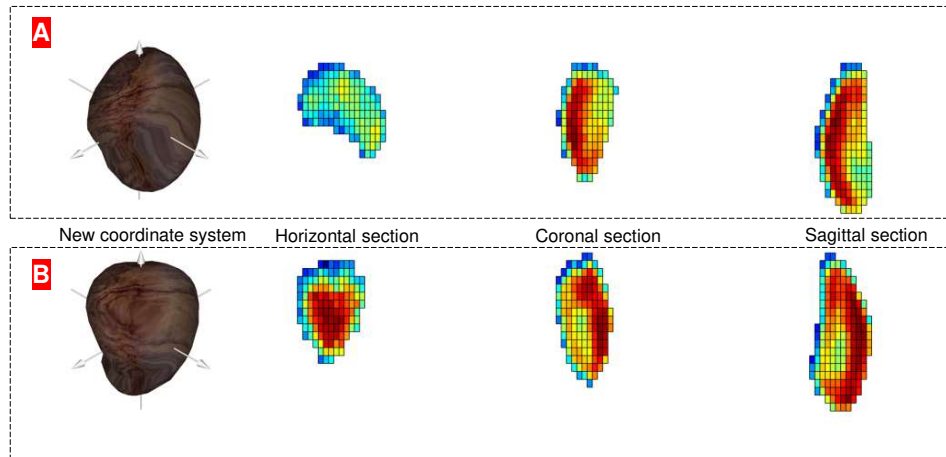


FIG. 3: Any section along any axis. (A) Sesamoid (left foot). (B) Sesamoid (right foot). The new coordinate system, based on the standardized body coordinate system, is the sesamoid posture after the anteroposterior axis rotated 12 degrees. Coronal section and sagittal section are the sections of new coordinate system with the same COM.

The COM and principal axis of inertia of the sesamoid are determined by the morphology and structure of the sesamoid. The body coordinate system built upon the sesamoid's COM and axes of inertia is unique. Its morphology and structure are irrespective of the posture and axis when the sesamoid is scanned. Therefore, the unique body coordinate can quantitatively describe the CSI of sesamoid (any cross section along any axis). See *Movie S2*.

- 
- [1] <http://www.mitk.net/index.html>
  - [2] M.J. Ciarelli, S. A. Goldstein, J. L. Kuhn, D. D. Cody, M. B. Brown, J. Orthop. R. 9, 674 (1991).
  - [3] J. Y. Rho, M. C. Hobatho, R. B. Ashman, Med. Eng. Phys. 17, 347(1995).
  - [4] Y. F. Fan, Y. B. Fan, Z. Y. Li, C. S. Lv, B. Zhang, PloS One 7, e32926 (2012).