

Fabrication of fully automated setup for high temperature thermal conductivity measurement

Ashutosh Patel* and Sudhir K. Pandey†

School of Engineering, Indian Institute of Technology Mandi, Kamand 175005, Himachal Pradesh, India

(Dated: June 11, 2022)

In this work, we report the fabrication of fully automated experimental setup for high temperature thermal conductivity (κ) measurement. Steady state based axial heat flow technique is used for κ measurement. Heat loss is measured using parallel thermal conductance technique. Simple design, light weight and small size sample holder is developed by using a thin heater and limited components. Low heat loss value is achieved by using small cross section very low thermal conductive insulator block. Power delivered to the heater is measured accurately by using 4-wire technique and for this, heater is developed with 4-wire. Program is built on graphical programming language "LabVIEW" to automate the whole measurement process. This setup is validated by using $Bi_{0.36}Sb_{1.45}Te_3$, polycrystalline bismuth, gadolinium, and alumina samples. The data obtained for these samples are found to be in good agreement with the reported data and the maximum deviation of 6 % in the value κ are observed. This maximum deviation is observed with the gadolinium sample. We also report the thermal conductivity of polycrystalline tellurium from 320 K to 550 K and the non monotonous behavior of κ with temperature is observed.

I. INTRODUCTION

In current scenario, finding new sources for power generation are very important to fulfill the demand of electricity. Various renewable energy sources like solar PV cell, wind turbine, etc. are discovered and currently are using at large scale. Various other sources have also been discovered but are in development stage. Power generation through thermoelectric (TE) materials is one of them. Research is going on in this field to find good TE materials, specially for high temperature application. A good TE material should have a high Seebeck coefficient (α), low thermal conductivity (κ), and high electrical conductivity (σ). Combining all three parameters, a common term introduced called figure of merit defined as $ZT = \alpha^2 \sigma T / \kappa$, where T is the mean temperature across TE material. A good TE material should have high ZT value. All three parameters have to be measured to find the value of ZT . Among these three properties, κ is most difficult to measure experimentally.

Various available methods to measure κ are discussed in Ref. [1]. Measurement methods are categorized in two parts, steady-state based and non-steady state based methods. Non steady state based methods like laser pulse method, 3ω -method, etc. require very less time for measurement as compared to steady state based methods. Laser pulse method is used most widely in commercial instruments but it is limited to thin disc sample only.² 3ω -method is also a popular technique, where time delay between the heating of the sample and the temperature response is used to measure κ .³ Steady state based methods, like Axial heat flow technique, comparative technique, guarded heat flow technique, heat-flow meter technique, etc. are also used for measurement.¹ The benefit of steady state based methods are that, setup cost is low and samples with various shapes and dimensions can be characterized. Steady state based axial heat flow method is a basic technique in which one-dimensional Fourier's

law of thermal conduction equation is used to measure κ , which can be written as

$$\kappa = \frac{\dot{Q}_s}{A} \cdot \frac{\Delta T}{l}$$

Where \dot{Q}_s is the net heat flow through the sample. A , l , and ΔT are cross-sectional area of the sample, its thickness, and temperature difference across it, respectively. The main difficulty with this method is to measure the heat flow through the sample (\dot{Q}_s) accurately. This difficulty arises due to the undefined amount of heat loss during the heat transfer process by conduction, convection and radiation, which becomes important at high temperature.⁴ Various instruments are reported based on this method for the measurement of κ . *Muto et al.*⁵ used a heat flux sensor to measure \dot{Q}_s . *Amatya et al.*⁴ used comparative axial heat flow technique to measure \dot{Q}_s . In this technique, a sample is sandwiched between standard reference materials of known κ . Radiation losses are taken care by using some physical model. *Zawilski et al.*⁶ discussed a different approach to find the value of \dot{Q}_s , called parallel thermal conductance technique. In this technique, heat loss is measured by running the instrument without sample. This heat loss data is act as a baseline for κ measurement. Using this technique, they performed measurement from 12 K to room temperature. Further, *Dasgupta et al.*⁷ used this technique and performed measurement at high temperature. They developed a model in which by using heat loss data the value of heater constant was obtained. The value of heater constant depends on the specific heat, mass, etc. of the heater block. Further this heater constant is used for thermal conductivity measurement. This technique was implemented by using a lifting mechanism for the heater. Heater heats the side wall of copper rod and this copper block heat the sample. They put the sample holder in a glass chamber and kept it in a furnace.

In the above discussed instruments, various techniques

are used to measure \dot{Q}_s . Use of heat flux sensor limit the measurement temperature, as it is not suitable for high temperature applications. Similarly, comparative axial heat flow technique requires standard reference materials and number of temperature sensors. Radiation loss at high temperature is estimated by using some physical model. To get a more accurate result by using the method discussed by *Zawilski et al.* heat loss should be as low as possible. Use of bulk heater causes more heat loss as it heat side walls and also high temperature surface is directly exposed to chamber environment.

In this work, we have developed a fully automated low cost setup for high temperature κ measurement. The Parallel thermal conductance technique is used for \dot{Q}_s measurement. We implemented this technique in such a way that, there is no requirement of furnace and lifting mechanism for heater. Heat loss is minimized by using low thermal conductive gypsum as an insulator block with small cross sectional area. Thin heater is built to heat the sample. Its negligible surfaces are exposed which also minimized the heat loss. The use of thin heater and small cross section insulator block minimizes the size of the sample holder assembly. This small size sample holder requires limited components, which make it lightweight. Thin heater is built with 4-probe contacts, by which we can measure accurate power delivered to the heater. LabVIEW based program with computer interfacing system is used to control the whole measurement process. $Bi_{0.36}Sb_{1.45}Te_3$, polycrystalline bismuth, gadolinium and alumina samples are used to validate the instrument. The data collected on these samples are found to be in good agreement with the reported data. We also report the κ of the polycrystalline tellurium.

II. MEASUREMENT METHODOLOGY

For the measurement of κ , one-dimensional Fourier's law of thermal conduction equation is used. As discussed in previous section, accuracy of the measurement mainly depends on the value of \dot{Q}_s . For this we have to find the heat loss value, which includes heat loss from the side walls of the sample, hot side copper block surface and insulator block surfaces. Heat loss occurred from the side walls of the sample by means of conduction, convection and radiation to the air. At high vacuum ($\sim 10^{-5} mbar$), this heat loss is less than 1 % of the total heat flow through the sample, even at high hot side temperature (above $500^\circ C$), which can be ignored.⁴ Heat loss in this system is mainly due to the heat flow by conduction through the insulator block, which is used to support the heater. Radiation losses from the hot side copper block surface and insulator surface are also considerable and cannot be ignored, specially at high temperature. So, mainly heat loss is due to the insulator block by conduction and radiation and hot side copper block by radiation.

For our system, heat generated by the heater is mainly divided into two paths. One is through sample and an-

other is through the insulator block. Heat loss through insulator block is obtained by running the instrument without sample. This method to find heat loss is known as parallel thermal conductance technique.⁶ Hot side temperature with respect to vacuum chamber temperature (T_{HR}) is recorded along with the power supplied to the heater. Room temperature is maintained at 300 K throughout the measurement process. As power supply changes, T_{HR} changes. The heat loss measurement is performed till the temperature for which κ measurement has to be performed. A schematic diagram for heat flow without sample and with sample are shown in Fig. 1. During heat loss measurement, the heat generated in the heater flows through the insulator block by conduction (\dot{Q}_{ins}) and through the surface of insulator block and hot side copper block by radiation (\dot{Q}_{rad}). The heat loss equation can be written as

$$\dot{Q}_l = \dot{Q}_{ins} + \dot{Q}_{rad}$$

Where \dot{Q}_l is the net power delivered to the heater during heat loss measurement.

During κ measurement, the sample is placed in the sample holder so the heat generated by the heater flows through the insulator as well as the sample also. The heat flow equation with a sample can be written as

$$\dot{Q}_s \approx \dot{Q}_{in} - \dot{Q}_l$$

Where \dot{Q}_{in} is the net power delivered to the heater during κ measurement. This \dot{Q}_s obtained using above equation is used to calculate the κ of the sample.

III. MEASUREMENT SETUP

The actual photograph of the sample holder assembly is shown in Fig. 2, where different components are represented by numbers. The sample, **1**, is sandwiched between two copper blocks, **2** & **3**, of $8mm \times 8mm$ cross section and 2 mm thickness. The two K-type PTFE shielded thermocouples, **4** & **5**, of 36 swg are embedded at the center of each copper block. Measured ΔT includes the temperature gradient generated due to the 1 mm of each copper block, the thermal resistance of interface surfaces and the sample. This instrument is used to measure sample of low thermal conductivity (up to 30 W/m.K), where the temperature gradient due to copper block is negligible as it has a very high thermal conductivity (~ 400 W/m.K at room temperature). GaSn liquid metal eutectic is used at the interface surfaces. Its thin film is generated at the interface surfaces, which minimized the interface resistance. Its boiling temperature is ~ 1800 K, which make it suitable at very high temperature also. 13Ω thin heater, **6**, is used to heat the sample. It is made by winding 40 swg kanthal wire over the alumina sheet of $7mm \times 7mm$ cross section and 0.6 mm of thickness. Each end of kanthal wire is welded with two nickel wires of 35 swg by using gas welding. At

each end, one wire is used to supply power and another is used to measure the voltage across kanthal wire. This 4-wire technique measures exact power delivered to the heater. This heater is further coated with high temperature cement and fixed over gypsum insulator block, **7**, by using high temperature cement. Gypsum is used for this purpose as it has a very low thermal conductivity (~ 0.017 W/m.K at room temperature) and can be used at high temperature also.⁹ As heat loss is very important and should be as low as possible, we take gypsum insulator block of $8\text{mm} \times 8\text{mm}$ cross section area and 25 mm thickness. Further, we minimized its effective cross sectional area by making a cylindrical hole across its thickness. This insulator block is supported by a brass plate, **8**, by using high temperature cement. This brass plate is fixed with a brass rod, **9**, using screw, **10**. Insulator block, **11**, is also used at the cold side to insulate sample electrically and optimize ΔT across it. SS rod, **12**, of 6 mm diameter and 60 mm length are used to apply pressure over the sample. This rod is having round tip at its end, which allows self aligning of the surfaces at the interface of sample and copper blocks. This self aligning ensures uniform pressure throughout the sample cross section. Another brass plate, **13**, is used to screw this SS rod. It is fixed over a plane ss rod, **14**, of 6 mm dia. Use of two different rods provide two independent paths for heat flow. Brass rod is for heat loss and ss rod is for heat flow through the sample. Brass rod and plane SS rod are fixed on a SS flange. On this SS flange, hermetically sealed electrical connector is fixed to make electrical connections. PT-100 RTD is connected to this connector to measure the temperature of the chamber. This sample holder assembly is placed inside the vacuum chamber, which is made by using seamless SS pipe of 10 cm diameter and 30 cm in height. KF25 port is provided over the vacuum chamber at the bottom. This port is used to connect this chamber with vacuum pump. Diffusion vacuum pump is used to create a vacuum inside the chamber upto a level of $\sim 3 \times 10^{-5}$ mbar. At this level of vacuum, the conduction and convection loss through the air can be ignored because of this pressure thermal conductivity of air is 3×10^{-6} W/m.K. It is four order of magnitude lower than thermal conductivity at atmospheric pressure.⁸

Keithley 2604B (SMU) is used to supply power to the heater and Keithley 2002 multimeter (DMM) with 10 channel scanner card is used to measure various signals. SMU is used in constant current mode and based on the resistance of the heater, power value is defined. The complete measurement process is controlled by a LabVIEW program. Before starting the program, the user needs to fill sample details, control parameters (step increment in the power supply, number of data average, and T_H limit). Heat loss data is also inserted into the string. A fixed power is supplied at each data point and based on this power, a temperature is achieved. Steady state is ensured by monitoring T_H , continuously and its rate of change is obtained. A very small rate of change in T_H

($0.1^\circ\text{C}/\text{min}$) is defined as steady state criteria. Once the rate of change of T_H comes below this set value, measurement starts. The random error generated in the measurement is minimized by taking an average of multiple data. At each data point T_H , T_C voltage across the heater and current through it are measured. The value of \dot{Q}_l is obtained by interpolating heat loss data and then the value of \dot{Q}_s is obtained by subtracting \dot{Q}_l from heater power. At the next step, κ is calculated using all raw data. This κ is plotted online with \bar{T} and all raw data along with κ is exported to .xls file at a user defined directory. After completion of this measurement loop, SMU increases the power supply to heater by the value defined in control parameter. This measurement loop continues until T_H reaches a set value defined in control parameter. Once T_H reaches set value, SMU stops power supply to the heater.

IV. RESULTS AND DISCUSSIONS

A. Heat loss measurement

Heat loss values are very important for κ measurement. At the start of the measurement, the value of \dot{Q}_s is very low. A small error in its value can change the result by considerable magnitude. For better result, \dot{Q}_l should be as low as possible. It is achieved by using low thermal conductive gypsum as insulator block. Again, the heat loss value is further minimized by reducing the effective surface area. A plot of heat loss with respect to T_{HR} is shown in Fig. 3. This heat loss data follow cubic degree polynomial, which is obtained by fitting of heat loss data with respect to T_{HR} . Initially, ~ 3 mW/K of heat loss is observed. This rate is increasing with increase in temperature, which is due to the increase in thermal conductivity of gypsum and radiation losses.

B. Thermal conductivity measurement

We calibrated our instrument by using $\text{Bi}_{0.36}\text{Sb}_{1.45}\text{Te}_3$, polycrystalline bismuth, gadolinium, and alumina samples. Samples with various κ values (varies from 1.4 W/m.K to 25 W/m.K) have been used to show the flexibility of the instrument. $\text{Bi}_{0.36}\text{Sb}_{1.45}\text{Te}_3$ sample is extracted from commercially available thermoelectric generator (TEC1-12706). The composition of the sample is obtained by performing EDX analysis. We also report the thermal conductivity of polycrystalline tellurium, which is not available to the best of our knowledge.

To verify the reproducibility of data on this instrument, measurement performed two times successively on the same sample. The measured values of κ are averaged. This average value of κ along with error bar are shown for all samples.

The dimensions of the samples used for measurement are given in the table.

Sample name	Cross section type	Cross sectional area (mm^2)	Thickness (mm)
$\text{Bi}_{0.36}\text{Sb}_{1.45}\text{Te}_3$	Rectangular	1.96	1.6
Poly. Bismuth	Rectangular	63	12
Gadolinium	Rectangular	2.56	2.4
Alumina	Rectangular	24	20
Poly. Tellurium	Rectangular	53	10

Measured and reported values of κ for $\text{Bi}_{0.36}\text{Sb}_{1.45}\text{Te}_3$ sample, polycrystalline bismuth, gadolinium, and alumina samples at various temperatures are shown in Fig. 4 (a), 4 (b), 4 (c), and 4 (d), respectively. Variation in ΔT with \bar{T} are shown in their insets. Measurements on $\text{Bi}_{0.36}\text{Sb}_{1.45}\text{Te}_3$, polycrystalline bismuth, gadolinium, and alumina samples are performed upto the mean temperature of 502 K, 440 K, 590 K, and 578 K, respectively and at these temperatures, the value of ΔT are 90 K, 26 K, 42 K, and 120 K, respectively.

Initially, the value of κ for $\text{Bi}_{0.36}\text{Sb}_{1.45}\text{Te}_3$ sample is 1.57 W/m.K. It decreases with increase in temperature with very small rate till $\bar{T}=350$ K and reaches to 1.565 W/m.K. After this temperature, κ increases with temperature and reaches to 2.67 W/m.K at $\bar{T}=502$ K. *Ma et al.*¹⁰ reported κ , in which they used commercial BiSbTe ingot as a sample. They obtained composition of sample by comparing XRD data of the sample with available XRD data of $\text{Bi}_{0.5}\text{Sb}_{1.5}\text{Te}_3$. Our data show similar behavior compared to the reported data. From 320 K to 375 K, almost constant deviation of 0.22 W/m.K is observed. It increases with increase in temperature and reach to 0.5 W/m.K at 500 K. From the experimental data, we observed that by increasing the concentration of bismuth, the value of κ increases.¹¹ The concentration of bismuth in our sample is more compared to the sample used in the reference, which may be the reason that our κ value is higher than the reported value.

In Bismuth, the electron and phonon components of κ are comparable in magnitude above 120 K. The value of κ decreases with temperature due to decrease in the magnitude of electron and phonon components with temperature.⁷ For a polycrystalline bismuth sample, we observed similar behavior. The value of κ at $\bar{T}=318$ K is 7.57 W/m.K. It decreases with increase in temperature and reaches to 6.73 W/m.K at $\bar{T}=440$ K. Our data shows similar behavior with the data available in handbook.¹² A deviation of 0.3 W/m.K is observed at $\bar{T}=318$ K. This deviation decreases with temperature and reach to 0.08 W/m.K. The maximum deviation of 0.3 W/m.K is observed at 318 K, which is 4 % of the value of the κ .

Gadolinium is a rare earth metal, which shows increasing behavior in thermal conductivity with temperature after ferromagnetic-antiferromagnetic transition. Its ferromagnetic-antiferromagnetic transition occurs near room temperature.¹³ Our data of gadolinium sample measured from 315 K to 590 K shows increasing behavior. Its thermal conductivity at $\bar{T}=315$ K is 7.76

W/m.K, with temperature κ increases and at $\bar{T}=590$ K, its value is 10.9 W/m.K. Data of gadolinium sample are compared with the reported data available in Ref. [14]. In this reference, values of κ are obtained by using thermal diffusivity data. A deviation of 0.37 W/m.K is observed at 315 K. This deviation decreases with temperature and match closely with reference data at 340 K. At 475 K and 575 K, deviation of ~ 0.18 W/m.K is observed. The maximum deviation of 0.37 W/m.K is observed at 315 K, which is the 5 % of the reported value of the κ .

It is observed that the thermal conductivity of oxide material decreases with temperature.¹⁵ Our data for alumina show similar behavior. At $\bar{T}=310$ K, the value of κ is 27.1 W/m.K. It decreases with increase in temperature and reaches to 19 W/m.K at $\bar{T}=578$ K. For alumina sample, the values of κ from various sources are compared in Ref. [16]. Near room temperature its value lies in the range of 24 W/m.K and 30 W/m.K. It decreases with temperature and at 600 K its value lies between 13.5 W/m.K to 17.5 W/m.K. Our data show similar behavior compared to the reported data. At 310 K, our data lies in the range of the reported data and at 578 K, a deviation of 0.5 W/m.K is observed.

We also performed measurement over polycrystalline tellurium sample. It is taken from commercially available ingot. The values of κ for this sample with respect to \bar{T} are shown in Fig. 5. Variation in ΔT along with \bar{T} is shown in the inset of Fig. 5. At $\bar{T}=320$ K, ΔT is 3 K. ΔT increases and reaches to 90 K at $\bar{T}=550$ K. At $\bar{T}=320$ K, κ is 3.14 W/m.K. It increases with increase in temperature with a very small rate till $\bar{T}=400$ K, where its value is 3.3 W/m.K. After this temperature, its value decreases and reached to 2.95 W/m.K at $\bar{T}=550$ K.

Thermoelectric properties of polycrystalline tellurium with different charge carrier concentrations are studied and measured with by *Lin et al.*¹⁷. Sample is prepared by melting and followed by quenching in cold water and annealing. charge carrier concentration tuned by using dopants including phosphorus (P), arsenic (As), antimony (Sb) and bismuth (Bi). With different carrier concentration, the value of κ changes. Its maximum value near room temperature is ~ 2.25 W/m.K for highest carrier concentration and with temperature it decreases.

A wide variety of samples are used for measurement, whose κ value lies in the range of 1.5 W/m.K to 27 W/m.K. At high temperature, previously only *Dasgupta et al.*⁷ used parallel conductance technique for κ measurement. They performed high temperature measurement with samples only upto a value of 2.5 W/m.K. We used this technique and performed measurement upto a value of 27 W/m.K. We also found reproducibility of our data with a maximum error bar of ± 3 %. We also observed a maximum deviation of 6 % compared to the reported data.

V. CONCLUSION

In this work, we have developed simple, low cost and user friendly setup for high temperature thermal conductivity measurement. Good reproducibility of measured data is observed with a maximum error bar of $\pm 3\%$. The parallel thermal conductance technique is used for measurement and data measured with better accuracy than previously reported literature using this method. We measured κ at high temperature using samples upto a value of 27 W/m.K. This is achieved by low heat loss and accurate power measurement by using 4-wire technique. Thin heater, simple design and limited components, make this setup light weight and small in size. The LabVIEW based program makes the whole measure-

ment process fully automated. This setup is validated by using $Bi_{0.36}Sb_{1.45}Te_3$, polycrystalline bismuth, gadolinium, and alumina samples. The measured data were found in good agreement with the reported data and observed a maximum deviation of 6 % with the gadolinium sample, which indicate that this instrument is capable to measure κ with fairly good accuracy. We also reported the thermal conductivity of polycrystalline tellurium and observed the non monotonous behavior of κ with temperature.

VI. ACKNOWLEDGEMENTS

The authors acknowledge R S Raghav and other workshop staffs for their support in the fabrication process of vacuum chamber and sample holder parts.

-
- * Electronic mail:ashutosh_patel@students.iitmandi.ac.in
 † Electronic mail:sudhir@iitmandi.ac.in
- ¹ X. C. Tong, "Characterization Methodologies of Thermal Management Materials," in *Advanced Materials for Thermal Management of Electronic Packaging*. **vol. 30**: Springer New York, pp. 59-129 (2011).
 - ² https://en.wikipedia.org/wiki/Thermal_conductivity_-_measurement
 - ³ H. Wang, M. Sen, Int. J. Heat Mass Transfer **52**, 2102-2109 (2009).
 - ⁴ R. Amatya, P. M. Mayer, and R. J. Ram, Rev. Sci. Instrum. **83**, 075117 (2012).
 - ⁵ A. Muto, D. Kraemer, Q. Hao, Z. F. Ren and G. Chen, Rev. Sci. Instrum. **80**, 093901 (2009).
 - ⁶ B. M. Zawilski, R. T. L. IV, and T. M. Tritt, Rev. Sci. Instrum. **72**, 1770 (2001).
 - ⁷ T. Dasgupta, A. M. Umarji, Rev. Sci. Instrum. **76**, 094901 (2005).
 - ⁸ A. Trowbridge, Phys. Rev. **2**, 58 (1913).
 - ⁹ S. L. Manzello, S. Park, T. Mizukami, and D. P. Bentz, Measurement of thermal properties of gypsum board at elevated temperatures., *Proceedings of the Fifth International Conference on Structures in Fire (SiF08)*, pp. 656-665 (2008).
 - ¹⁰ Y. Ma, Q. Hao, B. Poudel, Y. C. Lan, B. Yu, D. Z. Wang, G. Chen, and Z. F. Ren, Nano Lett. **8**, 2580 (2008).
 - ¹¹ J. Jianga, L. Chena, S. Baia, Q. Yaoa, Q. Wang, Journal of Crystal Growth, **277**, 258263 (2005).
 - ¹² C. Y. Ho, R. W. Powell and P. E. Liley, J. Phys. Chem. Ref. Data **3**, 1 (1974).
 - ¹³ W. J. Nellis, "Thermal conductivity and Lorenz function of gadolinium, terbium, and holmium single crystals" (1968). Retrospective Theses and Dissertations. Paper 4618.
 - ¹⁴ C. Meist, A. K. Froment and D. Moulinier, J. Phys. D: Appl. Phys **26**, 560-2 (1993).
 - ¹⁵ W. D. Kingery, J. Francl, R. L. Coble, and T. Vasilos, J. Am. Ceram. Soc., **37**, 107-10 (1954).
 - ¹⁶ P. Auerkari, Mechanical and physical properties of engineering alumina ceramics, Technical Research Center of Finland (VTT), Finland (1996).

- ¹⁷ S. Lin, W. Li, Z. Chen, J. Shen, B. Ge, and Y. Pei, Nat. Commun., **7**, 10287, (2016).

VII. FIGURE CAPTIONS

FIG. 1. Schematic diagram of heat flow without sample and with sample.

FIG. 2. Photograph of the sample holder.

FIG. 3. Heat loss values at different T_{HR} .

FIG. 4. (a) Thermal conductivity of $Bi_{0.36}Sb_{1.45}Te_3$ at different \bar{T} (variation in ΔT with \bar{T} is shown in the inset), (b) Thermal conductivity of polycrystalline bismuth at different \bar{T} (variation in ΔT with \bar{T} is shown in the inset), (c) Thermal conductivity of gadolinium at different \bar{T} (variation in ΔT with \bar{T} is shown in the inset), (d) Thermal conductivity of alumina at different \bar{T} (variation in ΔT with \bar{T} is shown in the inset).

FIG. 5. Thermal conductivity of polycrystalline tellurium at different \bar{T} (variation in ΔT with \bar{T} is shown in the inset).

FIG. 1.

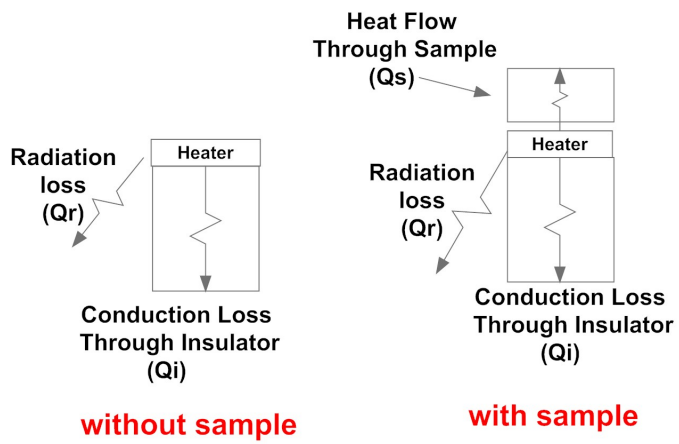
Heat Flow Diagram:

FIG. 2.

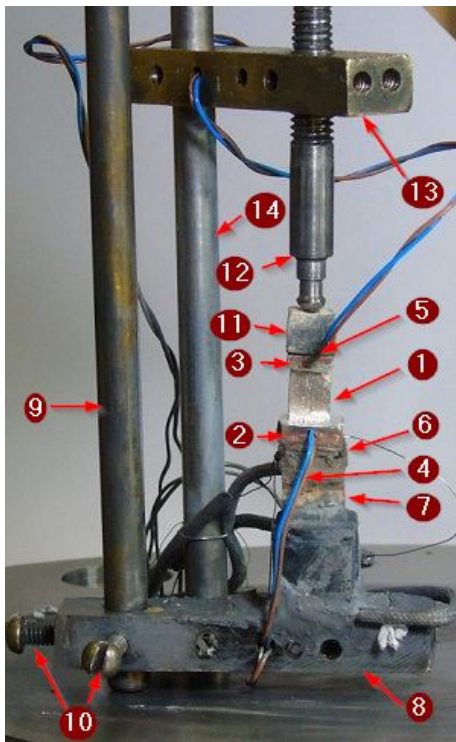


FIG. 3.

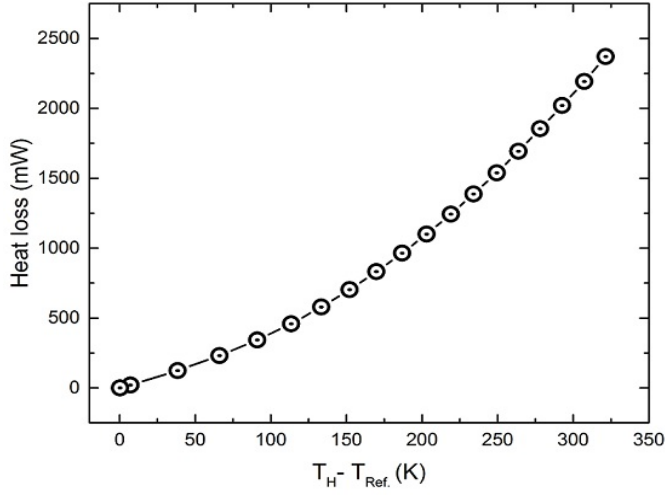


FIG. 4.




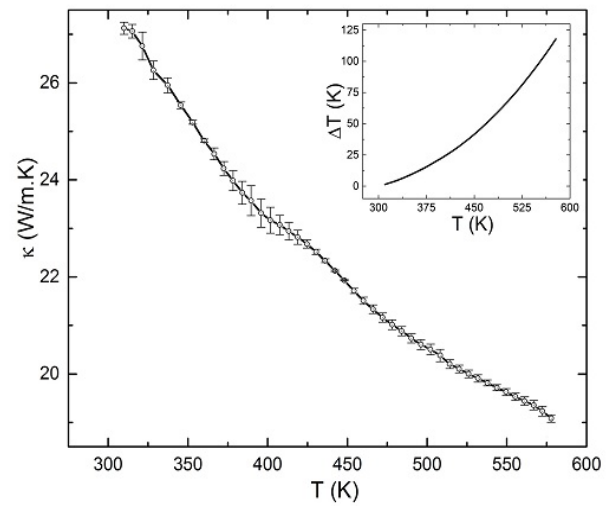


FIG. 5.

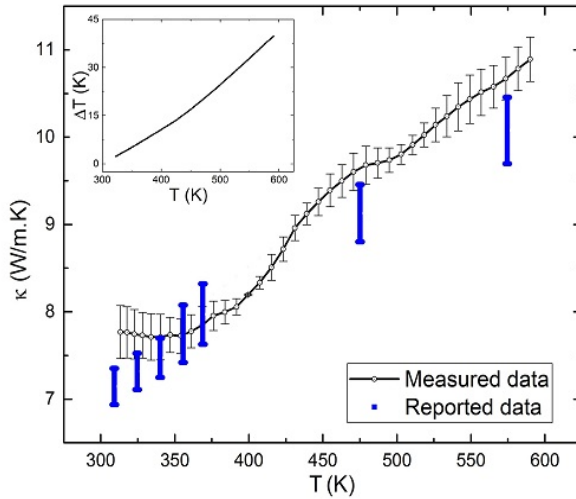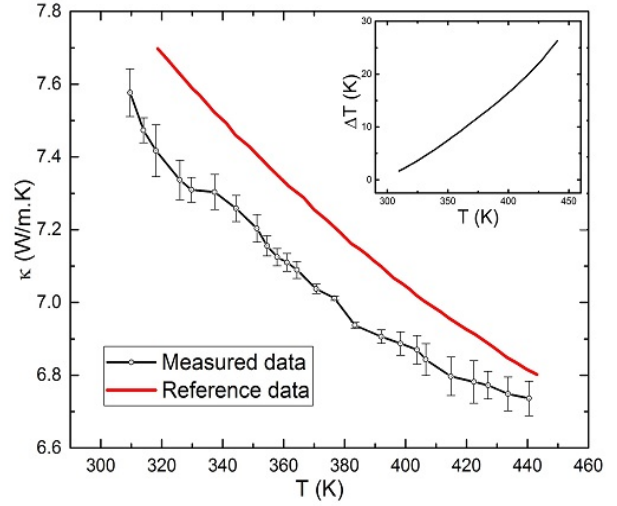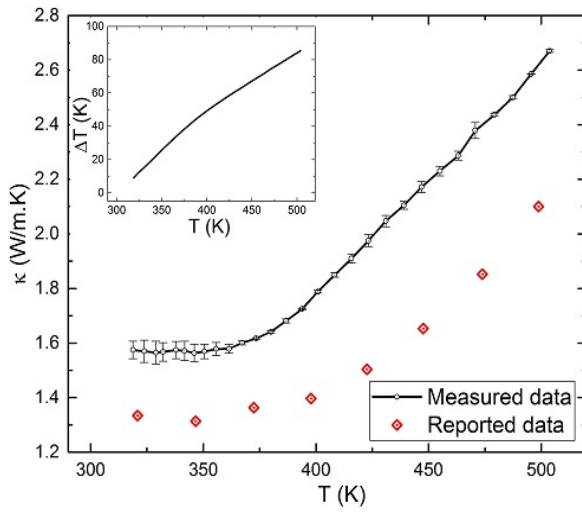