

Moore's Law in Single Cell Transcriptomics

Valentine Svensson^{1,2}, Roser Vento-Tormo², Sarah A Teichmann^{1,2}

Author Affiliations

¹European Molecular Biology Laboratory, European Bioinformatics Institute (EMBL-EBI), Hinxton, Cambridge, UK.

²Wellcome Trust Sanger Institute, Hinxton, Cambridge, UK.

Abstract

The ability to measure the transcriptome of single-cells has only been feasible for a few years, but is becoming an extremely popular assay. While many types of analysis and questions can be answered using single cell RNA-sequencing, of prime interest is the ability to investigate what cell types occur in nature. Unbiased and reproducible cataloging of distinct cell types require large numbers of cells to be sampled. Technological development and improvement of protocols has exponentially scaled the size of single cell RNA-seq studies, much faster than Moore's Law. In this perspective we will illustrate the steps that made this growth happen, and what that has meant for our ability to define cell types.

Introduction

Cells are the fundamental units of life, and cell differentiation allows the generation of complex multicellular organisms. This variety of cells states is made possible because cells exhibit different identities that are determined by their own nature (cell-intrinsically, through their development) and nurture (environment). The molecular mechanisms that regulate cell state and function are of fundamental interest, as they reflect underlying phenotype, as well as informing on developmental origin and tissue context.

Regulation of cell state and function occurs significantly at the mRNA level, with transcription factors defining a cell's transcriptional "program". In light of this, the abundances of different RNAs within cells are representative of functionally relevant cellular states. The ability to measure RNA in cells ties together cell biology, in terms of cellular phenotypes, with molecular biology, in terms of regulation of function.

Shortly after the advent of RNA sequencing as a biological tool¹, researchers started working on adapting the technology to single cell assays. Imaging-based assays such as in situ hybridization for RNA, or immunostaining for proteins had long revealed that population level averages were not representative of individual cell states².

Over the last few years, many sensitive and accurate single-cell RNA-sequencing (scRNA-Seq) protocols have been introduced³. In this Perspective piece, we show that single cell transcriptome experiments have grown exponentially, up to hundreds of thousands cells per study, in less than a decade (Figure 1). We also highlight what this has meant for the investigation of cell states, cell types and other sources of cell-to-cell variation.

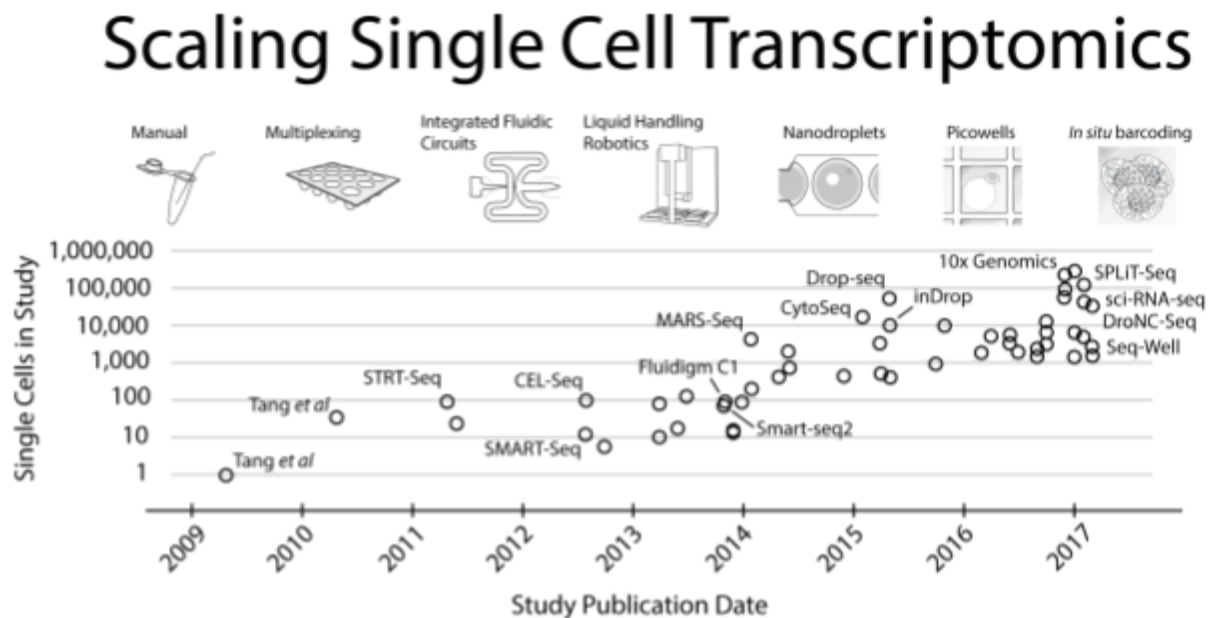


Figure 1: Scaling Single Cell Transcriptomics. Since the advent of single cell RNA sequencing, key technologies (top) have advanced the scale of single cell RNA-seq studies. Representative studies over the years are illustrated with their publication dates compared to their sample sizes (bottom). Key technological advances are named and annotated in the figure close to the corresponding study.

From tens to hundreds of cells

The increase in multiplexing of samples is driven by two factors: (i) reduction in the volume of reagents per sample and (ii) parallel processing of samples.

In the first single-cell RNA-seq study published in 2009 by Tang *et al*⁴, a single 4-cell stage blastomere was manually isolated using a glass capillary. The entire RNA-sequencing procedure was performed on the cell individually. The motivation was to make the most of precious embryonic samples compared to microarray techniques. The original method by Tang *et al* requires multiple PCR tubes for each cell, and a gel purification step^{5,6}.

Since it was known that intrinsic transcriptional variation might cause measurement issues, a technique called STRT-Seq was developed and presented in a study of Islam *et al*⁷ to multiplex up to 96 samples at once, in order to characterize the transcriptional landscape of MEFs and mESCs. At the time, the value of multiplexing strategies started being recognized as a means to allow higher throughput of sequencing experiments in general⁸. In the technique, the 96 cells are added to individual wells, where cell specific barcodes are added using the template-switching mechanism of reverse transcriptase during cDNA generation independently for each sample. The material from each well is then pooled together before it

is amplified by PCR. From the results of Islam *et al*⁷, it was clear that while there was a great degree of intrinsic variation between cells of the same type, embryonic stem cells (mESC) could be distinguished from MEF cells based on the transcriptome through clustering, illustrating the notion of a transcriptional cell state.

The ability to distinguish cell types by single cell transcriptomics was also demonstrated by Hashimshony *et al*⁹, where they introduced CEL-seq technique. As with STRT-seq, cells in individual wells are barcoded during cDNA synthesis and pooled. But rather than using PCR, the material is amplified by in vitro transcription (IVT).

Illumina, Inc. introduced highly multiplexed procedures for Illumina sequencing using dual-index barcoding in the form of the transposase-based 'tagmentation' with the Nextera XT kits¹⁰. This meant that multiplexing cells could be simplified by commercially available kits, allowing up to 192 samples per Illumina sequencing lane. The SMART-seq technique was introduced as a single-cell sequencing technology by Ramskold *et al*¹¹, giving full length coverage using the SMART template switching technology. The data generated was more familiar to users of traditional RNA-sequencing, making bioinformatics processing easier. This method made use of the Nextera kit for multiplexing and library generation. The Smart-seq technique made scRNA-seq experiments widely accessible by becoming readily available as the SMARTer kit from Clontech¹². A drawback compared to STRT- and CEL-seq is that amplification and library generation must be generated for all the cells individually, adding labour to the process.

While readily available and easy to use (Smart-seq has a smaller number of experimental steps than STRT-Seq or CEL-Seq), the cost per volume of the SMARTer kit was prohibitive to scale up the number of samples to the numbers in the STRT-seq and CEL-seq applications. Indeed, none of the published studies using this technique in microwells processed more than 20 cells^{11,13,14}. While it was noted that discrete cell states could be identified¹⁴, it was clear that more statistical power from more samples would be helpful. Two parallel strategies emerged to tackle this: reducing the reaction volume per cell, and adaptation of cheaper reagents for a similar goal.

The introduction of the Fluidigm C1 system improved the reagent cost, by letting the reaction occur in nanoliter chambers of an integrated fluidic circuit (IFC)¹⁵. The system also simplified cell isolation, as users can simply load a cell suspension into the system. Cells would be automatically captured in 96 chambers of the IFC, and processed material could later be transferred to microwell plates^{16,17}. Recently, a higher throughput version of the IFC allowing up to 800 cells to be captured was announced¹⁸. Reaching up to ~100 samples allowed researchers to estimate variance to and decompose it into biological and technical noise^{17,19}.

The Linnarsson group also modified their STRT-Seq protocol to be compatible with the C1 IFCs, which in addition reduced the reagent volume per cell, reduced the labour of manually isolating cells²⁰. This eventually enabled the group to survey and catalog neuronal subtypes in mouse cortex by investigating 3,005 cells from 67 individual mice²¹, after the addition of robotic automation for library preparation. The large number of cells from multiple mice allowed them to develop BackSPIN, a stable algorithm to cluster cells into neuronal cell types.

The group of Rickard Sandberg released the Smart-seq2 protocol, a refined version of Smart-seq which used less expensive, off the shelf reagents^{22,23}, which also allowed smaller volumes in individual wells. The microwell based Smart-seq2 protocol was more appropriate when cells could not be provided in a dense suspension, as is the case e.g. when studying early embryonic cells²⁴.

Researcher also noted that some cell types had troubles with the capture efficiency of the C1 IFC chips (for example, in the study by Zeisel *et al*²¹ the average successful cell capture rate was 41%). With the lower reagent price per volume, the Smart-seq2 protocol gave an alternative with greater control, though at the cost of labour. The strategy of relying on off-the-shelf reagents and not needing the expensive IFC reduced the price per cell. The dominant cost factor for Smart-seq2 at this point was the reliance on the Nextera kit for barcoding. The cost per cell was further reduced by in-house production of Tn5, a variant of the active enzyme in Illumina's Nextera XT fragmentation kit²⁵. The same Tn5 transposase was used by the Linnarsson lab to reduce the cost of the STRT-seq technology²⁰.

Surpassing thousands of cells

With well-calibrated flow cytometers, cells can quickly be isolated in individual wells of 96 or 384 microwell plates. Once cells are available as a single cell suspension, researchers can use this strategy to populate large numbers of plates with their cells of interest, the bottleneck becoming the processing of the plates. With this in mind, Jaitin *et al* modified the CEL-seq protocol to be compatible with robotic automation in MARS-seq²⁶. This allowed the Amit team to decrease the labour of processing plates filled with isolated cells, and scale up massively, investigating 4,000 cells in one study²⁶. Several laboratories since have set up automation procedures for some standard protocol of choice^{21,27}, and the refined CEL-seq2 method²⁸ was automated in form of the SORT-seq method²⁹. The large number of samples allowed the researchers to develop a probabilistic mixture model and assign cells to immune cell types without known markers²⁶.

In CEL-seq and derivative protocols^{9,26,28,30}, it had been demonstrated that as long as you had isolated and barcoded the cDNA material, the following steps of amplification and library preparation could be done in a single unit. The throughput bottleneck was pinpointed to two major factors: isolation of cells, and the ability to generate enough multiplexing barcodes to investigate large numbers of cells in parallel. Methods had emerged to randomly capture and manipulate individual cells in nanoliter droplet emulsions³¹. The challenge of creating cDNA and barcode the material in the individual droplets was solved by Klein *et al*³² and Macosko *et al*³³ by delivering beads coated with barcoded sequences into the droplets in the inDrop and Drop-seq methods, respectively.

The reverse droplets have miniscule volumes of about 1 nL per cell, further reducing the reagent cost per cell. However, to avoid capturing two cells in a single droplet, the random isolation must be rate limited. This means the efficiency of cells per unique bead is low (on the order of 5-10%^{32,33}), and the barcode space must be sufficiently large to allow for a great number of unused sequences. This is a challenge considering the molecular properties of the barcodes need to be accounted for during the barcode design³⁴.

A benefit of the droplet microfluidics is that it is simple to manufacture the components. The McCarroll lab has made plans and details for the Drop-seq system available online at <http://mccarrolllab.com/dropseq/>. This spreads the use of the technology, and has enabled researchers to customize it to their needs. Still, optimization and experience with microfluidics is needed for optimal results. The company 10x Genomics commercialised the material required (device and reagents) to perform the inDrop method, spreading the technology worldwide³⁵. In their implementation, up to 8 independent cell pools can be processed simultaneously, allowing the parallelisation of several experiments in a single run. The technology was demonstrated in a massive study of 250,000 cells³⁶. Recently, Illumina and Bio-Rad also announced nanodroplet based single cell system³⁷.

An alternative strategy for massively parallel cell isolation is to deposit beads into picoliter wells, and randomly load them with cells^{38–40}. Beside decreasing the reaction volume, the picowell systems are easier to control, and more portable. This allows more rapid collection of fresh cells in e.g. clinical settings, further increasing the ability to make single cell observations.

A related strategy by Vickovic *et al*⁴¹ consists of FACS sorting cells onto the surface of an array with attached barcoded Poly(dT) capture probes. The throughput is estimated to be 10,000 cells over two days. Since cells are not randomly isolated but actively placed, the efficiency per cell should be better depending on the flow sorter.

On the horizon: hundreds of thousands to millions of cells

Further scaling of random isolation methods, through nanodroplets or picowells to analyse cells in parallel is limited by rate-limiting to avoid doublets of cells, as well as the number of synthesized barcodes. Two new methods have been demonstrated to overcome these issues, by making use of sequential *in situ* barcoding^{42,43}.

In the *in situ* barcoding approaches, cells are permeabilized and fixed with formaldehyde⁴² or methanol⁴³. The fixed cells are then divided into small fractions of reactions where in the first round cDNA is generated from RNA, and all cells in a fraction get a unique barcode. In subsequent steps, cells are pooled together, and re-divided into fraction where the new fractions get another barcode appended to the cDNA. The low probability of cells going together into the sequential fractions means each cell will get a unique sequence of barcodes, while the cells have actually never been isolated. In the study of Cao *et al*, this technique was applied to *C. elegans* which gave the researchers about 50 copies of each known cell type in the worm, providing a single-cell atlas of a whole animal⁴³. Cell types could be identified in an unbiased manner by first creating a low-dimensional representation through, then applying the Density Peak clustering algorithm⁴⁴. The high copy number of the cell types allowed the low-dimensional representation to be well represented with distinct groups.

Currently, the major cost limitation is library preparation and cDNA sequencing. While it has been shown that to identify cell types and regulatory networks, relatively shallow sequencing

suffice^{45,46}, the sequencing cost is still prohibitive even at shallow depths when reaching hundreds of thousands of cells. Recent announcements promises slightly cheaper sequencing from higher throughput⁴⁷, some radical change in sequencing technology might be needed.

Beyond sequencing, the limitation is obtaining cells. In some cases, cells might not be easy to isolate, and in many interesting cases, it is difficult to get cells into a single cell suspension. In line with this, two laboratories in parallel adapted the single-cell transcriptomics to nuclei isolation, allowing the work with tissues where harsh dissociation protocols will damage RNA integrity^{48–50}. Moreover, stored material can be used in nuclear single-cell RNA analysis, as it is compatible with fixed cells. Recent work has shown that cells can be preserved prior to preparation for single cell RNA-sequencing, increasing the ability to gather usable material without reducing the complexity of the transcriptome^{51,52}.

A particularly illustrative example of the scaling which have happened over the year is from the Regev lab scaled 18 cells¹⁴ (microwell plates) to 2,000 cells⁵³ (microfluidics) to 200,000 cells⁵⁴ (nanodroplets) of the same cell type over the last 4 years. The field of transcriptomics have long been cursed by the so called “large P small N problem”, where the number of observations (cells) is much smaller than the number of variables (genes). With hundreds of thousands of cells, each expressing up to 10,000 genes, this is no longer the case, and in the coming years we will see a lot of interesting results from this.

Recently, single cell studies have started including artificial perturbations of the system, from which direct regulatory information can be gained using relatively simple linear models^{54–57}. It has been known that simple analysis methods with a lot of data perform better than complex methods with small amounts of data⁵⁸. A recent good example of this is in the work of Esteva *et al*, where a gold standard data set of 1,000 images was increased to one with 100,000, images, allowing the researchers to train a neural network which beat trained dermatologists in clinical classification of skin lesions⁵⁹. When studies reach millions of cells, even rare cell types will be identifiable without issue, and we will have an unbiased view on transcriptional diversity.

Author Contributions

V.S., R.V-T., and S.A.T. wrote the manuscript. V.S. prepared the figure.

Competing Financial Interests

The authors declare no competing financial interests.

Acknowledgements

The authors would like to thank Zheng-Shan Chong and Kedar Nath Natarajan for reading the manuscript and giving helpful feedback.

Figure Legends

Figure 1:

Scaling Single Cell Transcriptomics. Since the advent of single cell RNA sequencing, key technologies (top) have advanced the scale of single cell RNA-seq studies. Representative studies over the years are illustrated with their publication dates compared to their sample sizes (bottom). Key technological advances are named and annotated in the figure close to the corresponding study.

References

1. Mortazavi, A., Williams, B. A., McCue, K., Schaeffer, L. & Wold, B. Mapping and quantifying mammalian transcriptomes by RNA-Seq. *Nat. Methods* **5**, 621–628 (2008).
2. Raj, A. & van Oudenaarden, A. Nature, nurture, or chance: stochastic gene expression and its consequences. *Cell* **135**, 216–226 (2008).
3. Svensson, V. *et al.* Power analysis of single-cell RNA-sequencing experiments. *Nat. Methods* (2017). doi:10.1038/nmeth.4220
4. Tang, F. *et al.* mRNA-Seq whole-transcriptome analysis of a single cell. *Nat. Methods* **6**, 377–382 (2009).
5. Sasagawa, Y. *et al.* Quartz-Seq: a highly reproducible and sensitive single-cell RNA sequencing method, reveals non-genetic gene-expression heterogeneity. *Genome Biol.* **14**, R31 (2013).
6. Tang, F. *et al.* RNA-Seq analysis to capture the transcriptome landscape of a single cell. *Nat. Protoc.* **5**, 516–535 (2010).
7. Islam, S. *et al.* Characterization of the single-cell transcriptional landscape by highly multiplex RNA-seq. *Genome Res.* **21**, 1160–1167 (2011).
8. Kozarewa, I. & Turner, D. J. 96-plex molecular barcoding for the Illumina Genome Analyzer. *Methods Mol. Biol.* **733**, 279–298 (2011).
9. Hashimshony, T., Wagner, F., Sher, N. & Yanai, I. CEL-Seq: single-cell RNA-Seq by

- multiplexed linear amplification. *Cell Rep.* **2**, 666–673 (2012).
10. Illumina, Inc. Illumina Introduces Nextera® XT DNA Sample Preparation Kits. (2012).
 11. Ramsköld, D. *et al.* Full-length mRNA-Seq from single-cell levels of RNA and individual circulating tumor cells. *Nat. Biotechnol.* **30**, 777–782 (2012).
 12. Clontech Laboratories, Inc. Clontech Laboratories, Inc. Releases the SMARTer® Universal Low Input RNA Kit for Sequencing. (2013).
 13. Marinov, G. K. *et al.* From single-cell to cell-pool transcriptomes: stochasticity in gene expression and RNA splicing. *Genome Res.* **24**, 496–510 (2014).
 14. Shalek, A. K. *et al.* Single-cell transcriptomics reveals bimodality in expression and splicing in immune cells. *Nature* **498**, 236–240 (2013).
 15. Fluidigm Corporation. FLUIDIGM COMPLETES C1™ SINGLE-CELL AUTO PREP SYSTEM. (2013).
 16. Wu, A. R. *et al.* Quantitative assessment of single-cell RNA-sequencing methods. *Nat. Methods* **11**, 41–46 (2014).
 17. Brennecke, P. *et al.* Accounting for technical noise in single-cell RNA-seq experiments. *Nat. Methods* **10**, 1093–1095 (2013).
 18. Fluidigm Corporation. FLUIDIGM INTRODUCES A NEW MEDIUM-CELL HIGH-THROUGHPUT INTEGRATED FLUIDIC CIRCUIT. (2016).
 19. Kim, J. K. *et al.* Characterizing noise structure in single-cell RNA-seq distinguishes genuine from technical stochastic allelic expression. *Nat. Commun.* **6**, 8687 (2015).
 20. Islam, S. *et al.* Quantitative single-cell RNA-seq with unique molecular identifiers. *Nat. Methods* **11**, 163–166 (2014).
 21. Zeisel, A. *et al.* Cell types in the mouse cortex and hippocampus revealed by single-cell RNA-seq. *Science* **347**, 1138–1142 (2015).
 22. Picelli, S. *et al.* Smart-seq2 for sensitive full-length transcriptome profiling in single cells. *Nat. Methods* **10**, 1096–1098 (2013).

23. Picelli, S. *et al.* Full-length RNA-seq from single cells using Smart-seq2. *Nat. Protoc.* **9**, 171–181 (2014).
24. Deng, Q., Ramsköld, D., Reinius, B. & Sandberg, R. Single-cell RNA-seq reveals dynamic, random monoallelic gene expression in mammalian cells. *Science* **343**, 193–196 (2014).
25. Picelli, S. *et al.* Tn5 transposase and tagmentation procedures for massively scaled sequencing projects. *Genome Res.* **24**, 2033–2040 (2014).
26. Jaitin, D. A. *et al.* Massively parallel single-cell RNA-seq for marker-free decomposition of tissues into cell types. *Science* **343**, 776–779 (2014).
27. Soumillon, M., Cacchiarelli, D. & Semrau, S. Characterization of directed differentiation by high-throughput single-cell RNA-Seq. *BioRxiv* (2014).
28. Hashimshony, T. *et al.* CEL-Seq2: sensitive highly-multiplexed single-cell RNA-Seq. *Genome Biol.* **17**, 77 (2016).
29. Muraro, M. J. *et al.* A Single-Cell Transcriptome Atlas of the Human Pancreas. *Cell Syst* **3**, 385–394.e3 (2016).
30. Velten, L. *et al.* Single-cell polyadenylation site mapping reveals 3' isoform choice variability. *Mol. Syst. Biol.* **11**, 812 (2015).
31. Mazutis, L. *et al.* Single-cell analysis and sorting using droplet-based microfluidics. *Nat. Protoc.* **8**, 870–891 (2013).
32. Klein, A. M. *et al.* Droplet barcoding for single-cell transcriptomics applied to embryonic stem cells. *Cell* **161**, 1187–1201 (2015).
33. Macosko, E. Z. *et al.* Highly Parallel Genome-wide Expression Profiling of Individual Cells Using Nanoliter Droplets. *Cell* **161**, 1202–1214 (2015).
34. Costea, P. I., Lundeberg, J. & Akan, P. TagGD: fast and accurate software for DNA Tag generation and demultiplexing. *PLoS One* **8**, e57521 (2013).
35. 10x Genomics, Inc. 10x Genomics Announces Commercial Availability of Single Cell 3'

Solution at the AACR Annual Meeting 2016. (2016).

36. Zheng, G. X. Y. *et al.* Massively parallel digital transcriptional profiling of single cells. *Nat. Commun.* **8**, 14049 (2017).
37. Illumina, Inc. Illumina and Bio-Rad Launch Solution for Single-Cell Genomic Sequencing to Enable Robust Study of Complex Diseases. (2017).
38. Christina Fan, H., Fu, G. K. & Fodor, S. P. A. Combinatorial labeling of single cells for gene expression cytometry. *Science* **347**, 1258367 (2015).
39. Gierahn, T. M. *et al.* Seq-Well: portable, low-cost RNA sequencing of single cells at high throughput. *Nat. Methods* (2017). doi:10.1038/nmeth.4179
40. Bose, S. *et al.* Scalable microfluidics for single-cell RNA printing and sequencing. *Genome Biol.* **16**, 120 (2015).
41. Vickovic, S. *et al.* Massive and parallel expression profiling using microarrayed single-cell sequencing. *Nat. Commun.* **7**, 13182 (2016).
42. Rosenberg, A. B. *et al.* Scaling single cell transcriptomics through split pool barcoding. *bioRxiv* 105163 (2017). doi:10.1101/105163
43. Cao, J. *et al.* Comprehensive single cell transcriptional profiling of a multicellular organism by combinatorial indexing. *bioRxiv* 104844 (2017). doi:10.1101/104844
44. Rodriguez, A. & Laio, A. Machine learning. Clustering by fast search and find of density peaks. *Science* **344**, 1492–1496 (2014).
45. Heimberg, G., Bhatnagar, R., El-Samad, H. & Thomson, M. Low Dimensionality in Gene Expression Data Enables the Accurate Extraction of Transcriptional Programs from Shallow Sequencing. *Cell Syst* **2**, 239–250 (2016).
46. Pollen, A. A. *et al.* Low-coverage single-cell mRNA sequencing reveals cellular heterogeneity and activated signaling pathways in developing cerebral cortex. *Nat. Biotechnol.* **32**, 1053–1058 (2014).
47. Illumina, Inc. Illumina Introduces the NovaSeq Series—a New Architecture Designed to

Usher in the \$100 Genome. (2017).

48. Habib, N. *et al.* Div-Seq: Single-nucleus RNA-Seq reveals dynamics of rare adult newborn neurons. *Science* **353**, 925–928 (2016).
49. Habib, N. *et al.* DroNc-Seq: Deciphering cell types in human archived brain tissues by massively-parallel single nucleus RNA-seq. *bioRxiv* 115196 (2017). doi:10.1101/115196
50. Lake, B. B. *et al.* Neuronal subtypes and diversity revealed by single-nucleus RNA sequencing of the human brain. *Science* **352**, 1586–1590 (2016).
51. Guillaumet-Adkins, A. *et al.* Single-cell transcriptome conservation in cryopreserved cells and tissues. *Genome Biol.* **18**, 45 (2017).
52. Alles, J., Praktijnjo, S., Karaiskos, N. & Grosswendt, S. Cell fixation and preservation for droplet-based single-cell transcriptomics. *bioRxiv* (2017).
53. Shalek, A. K. *et al.* Single-cell RNA-seq reveals dynamic paracrine control of cellular variation. *Nature* **510**, 363–369 (2014).
54. Dixit, A. *et al.* Perturb-Seq: Dissecting Molecular Circuits with Scalable Single-Cell RNA Profiling of Pooled Genetic Screens. *Cell* **167**, 1853–1866.e17 (2016).
55. Jaitin, D. A. *et al.* Dissecting Immune Circuits by Linking CRISPR-Pooled Screens with Single-Cell RNA-Seq. *Cell* **167**, 1883–1896.e15 (2016).
56. Adamson, B. *et al.* A Multiplexed Single-Cell CRISPR Screening Platform Enables Systematic Dissection of the Unfolded Protein Response. *Cell* **167**, 1867–1882.e21 (2016).
57. Datlinger, P. *et al.* Pooled CRISPR screening with single-cell transcriptome readout. *Nat. Methods* **14**, 297–301 (2017).
58. Halevy, A., Norvig, P. & Pereira, F. The Unreasonable Effectiveness of Data. *IEEE Intell. Syst.* **24**, 8–12 (2009).
59. Esteva, A. *et al.* Dermatologist-level classification of skin cancer with deep neural networks. *Nature* **542**, 115–118 (2017).

Scaling Single Cell Transcriptomics

