

# *SiMRX* - A Simulation toolbox for *MRX*

Technical Report for *SiMRX* version 1.3

Janic Föcke<sup>1,2</sup>

July 29, 2022

## Abstract

*SiMRX* is a *MRX* simulation toolbox written in *MATLAB* for simulation of realistic 2D and 3D Magnetorelaxometry (*MRX*) setups, including coils, sensors and current patterns. *MRX* is a new modality that uses magnetic nanoparticles (*MNP*) as contrast agent and shows promising results in medical applications, e.g. cancer treatment. Its basic principles were outlined in [2], further elaborated in [6], transferred into a rigorous mathematical model and analyzed in [3].

*SiMRX* is available at <https://gitlab.com/simrx/simrx/>.

## Acknowledgements and Contributions

This toolbox was inspired by earlier work of Peter Hoemmen<sup>3</sup>, Paul Koenigsberger<sup>3</sup> and Daniel Baumgarten<sup>3,4</sup>. In particular for the simulation of the 3D case, we recycled code fragments originally created by Peter Hoemmen<sup>3</sup> and Paul Koenigsberger<sup>3</sup>. Moreover, the 3D simulation step has been validated using a measured 3D dataset provided by Maik Liebl<sup>5</sup>. The phantoms were created as part of a collaboration with Daniel Baumgarten<sup>3,4</sup> and Peter Schier<sup>4</sup>. As part of the ongoing collaboration with Daniel Baumgarten<sup>3,4</sup> and Peter Schier<sup>4</sup> this code may receive updates and new features in the future.

Maik Liebl has been supported by the German Science Foundation (DFG) within the priority program SPP1681 (WI 4230/1-3).

This work has been supported by the German Science Foundation (DFG) within the priority program CoSIP, project CoS-MRXI (BA 4858/2-1, BU 2327/4-1).

**keywords**— magnetorelaxometry imaging; magnetic nanoparticles; modeling; *MATLAB* toolbox; *SiMRX*

## 1 Introduction

Many new and experimental treatment methods in medical applications use magnetic nanoparticles as a contrast agent. These particles allow for multiple different approaches ([1, 5]), however for the named methods the exact knowledge about the particle distribution is crucial. Here, Magnetorelaxometry (*MRX*) can be used to determine the amount of particles in a region of interest as shown in [2, 6]. Based on this approach an imaging technique called Magnetorelaxometry Imaging (*MRXI*) has been proposed.

The *SiMRX* toolbox provides a set of tools to model and simulate such an *MRX* system. It is based on the mathematical model developed in [3].

<sup>1</sup>Institute for Analysis and Numerics, Westfälische Wilhelms-Universität Münster (WWU), DE

<sup>2</sup>Applied Mathematics, Friedrich-Alexander Universität Erlangen-Nürnberg (FAU), DE

<sup>3</sup>Institute of Biomedical Engineering and Informatics, Technische Universität Ilmenau (TU Ilmenau), DE

<sup>4</sup>Institute of Electrical and Biomedical Engineering, Private University of Health Sciences, Medical Informatics and Technology (UMIT), Hall in Tirol, AT

<sup>5</sup>Physikalisch-Technische Bundesanstalt (PTB), Berlin, DE

## 1.1 Model

This following section (including notation) is part of [3], which, for interested readers, provides an in depth look in the analysis of the following operator. The magnetic field in  $w \in \Omega$  induced by a coil  $\alpha = (\varphi_\alpha, I_\alpha)$  is given by

$$\mathbf{B}_\alpha^{\text{coil}}: \Omega \rightarrow \mathbb{R}^3, \quad w \mapsto \vartheta \int_0^{L_\alpha} \varphi'_\alpha(s) \times \left( \frac{w - \varphi_\alpha(s)}{|w - \varphi_\alpha(s)|^3} \right) ds, \quad (1)$$

where  $L_\alpha$  is the length of the coil. The magnetization of the magnetic nanoparticles (MNP) after the reorientation process in  $w$  is described by

$$\mathbf{m}_\alpha: \Omega \rightarrow \mathbb{R}^3, \quad w \mapsto \frac{1}{3} \mathbf{B}_\alpha^{\text{coil}}(w) c(w), \quad (2)$$

where  $c$  is the desired particle distribution. The particle induced magnetic response from particles in  $w$  measured by a sensor  $\sigma = (\sigma_x, \sigma_n)$  is modeled by

$$\mathbf{B}_\alpha^{\text{meas}}: \Omega \times \Sigma \rightarrow \mathbb{R} \\ (w, \sigma) \mapsto \sigma_n \cdot \left( \left( \frac{3(\sigma_x - w) \otimes (\sigma_x - w)}{|\sigma_x - w|^5} - \frac{\mathbb{I}}{|\sigma_x - w|^3} \right) \mathbf{B}_\alpha^{\text{coil}}(w) \right). \quad (3)$$

In the end we receive the following forward operator for a coil  $\alpha$ :

$$\mathbf{K}_\alpha: \mathcal{L}^2(\Omega) \rightarrow \mathcal{L}^2(\Sigma), \quad c \mapsto \left[ \sigma \mapsto \int_\Omega \mathbf{B}_\alpha^{\text{meas}}(w, \sigma) c(w) d^3 w \right]. \quad (4)$$

## 1.2 Discretization

First we consider the 3D case: Here the conductor coil  $\varphi_\alpha$  is approximated by a set list of segments, with starting points  $a_k$  and ending points  $b_k$  for the  $k$ -th segment respectively. Then the magnetic field in  $w \in \Omega$  is [4]:

$$\mathbf{B}_{\alpha,k}^{\text{coil}}: \Omega \rightarrow \mathbb{R}^3 \\ w \mapsto \vartheta \frac{|a_k - w| + |b_k - w|}{|a_k - w||b_k - w|} \frac{(a_k - w) \times (b_k - w)}{|a_k - w||b_k - w| + (a_k - w) \cdot (b_k - w)}. \quad (5)$$

The MNP response (see equation (3)) still holds in the discrete case.

In the 2D case a coil simplification leads to the usage of variants of (3) for both coil and dipole response (see [3, section 4.2]).

## 2 Structure

The *SiMRX* toolbox is a modular toolkit that provides tools for MRX experiment setup configuration and simulation as well as visualization of data. *SiMRX* is capable of simulating synthetic or real setups and datasets.

The processes required for simulation are handled in a sequence of modular functions, i.e. the creation of voxel grids, the calculation of magnetic fields and ultimately the data acquisition at the

sensors. See section 3.2 for a detailed description of the simulation step and relate to Figure 1 for an illustration of the simulation process.

For the simulation of an MRX experiment we introduce a distinction between a `setup` and `config` file. The `setup` file includes information about the number of dimensions, coil position and shape (and, if necessary, orientation), sensor position and orientation, as well as intervals that define the region of interest. The `config` file includes information about the resolution of the phantom within the region of interest, coil current patterns, as well as information what subset of coils and sensors are actively used. Both files datasets are forwarded to the simulation script, that return a matrix representation of the given setup and configuration. Detailed information on the creation of `setup` and `config` file can be found in section 3.1.

Furthermore *SiMRX* provides useful tools that visualize MRX setups and internal states (section 3.4), and includes a tool for the creation of phantoms (section 3.3).

**Please note:** The *SiMRX* toolbox used a unified coordinate system, where spatial information is stored in a three dimensional vector  $(x, y, z)$ . Any position is given in relation to a zeropoint  $(0, 0, 0)$  and, if not specified further, is provided in [m].

If *SiMRX* is used to implement a 2D setup, please read the paragraph about *Conventions for 2D configurations* in the introduction of section 3.1.

## 3 Features

The *SiMRX* toolbox is separated in the following submodules (each in its own subfolder):

```
./configuration
./simulation
./phantom
./visualization
```

We give a short overview of each module in the following sections. For all provided functions, in detail information of syntax and features are available in each file header. This documentation is also available using the MATLAB help function (F1).

### 3.1 Configuration

This module provides a set of functions to create `setup` and `config` files.

**Conventions for 2D configurations:** *SiMRX* supports 2D and 3D systems. However, at its core, the toolbox is designed for 3D simulations. Therefore 2D setups are implemented as 3D setups with the following variances:

- all entities are on the same z-layer ( $z=0$  recommended)
- the region of interest in z direction is set accordingly (`setup.roi.z = [0,0]` recommended)
- the voxel resolution `config.res` in z direction is set to 1, e.g. `config.res = [x,y,1]`;
- coils *must be* represented as dipole magnets to be in line with section 1.2. This means that `coils.Segments` will be ignored, regardless of availability.

Also consider the provided 2D and 3D setups as part of *SiMRX* (see section 4.2.1).

---

**Listing 1: Create a new 2D setup**

---

```
1  %% create a new 2D setup
2  setup = [];
3
4  setup.dim = 2;
5
6  setup.roi.x = [-.05, .05];
7  setup.roi.y = [-.05, .05];
8  setup.roi.z = [ 0, 0];
9
10 coils = [];
11 coils = [coils, createEntityArray([ .055,.045,0], [ .055,-.045,0], [-1,0,0], [1,10,1])];
12 coils = [coils, createEntityArray([-0.055,.045,0], [-.055,-.045,0], [ 1,0,0], [1,10,1])];
13 coils = [coils, createEntityArray([.045, .055,0], [-.045, .055,0], [0,-1,0], [10,1,1])];
14 coils = [coils, createEntityArray([.045,-.055,0], [-.045,-.055,0], [0, 1,0], [10,1,1])];
15 setup.coils = coils;
16
17 sensors = [];
18 sensors = [sensors, createEntityArray([.050,.060,0], [-.050,.060,0], [-1,-1,0], 10)];
19 sensors = [sensors, createEntityArray([.050,.060,0], [-.050,.060,0], [ 1,-1,0], 10)];
20 setup.sensors = sensors;
21
22 visualizeSetup(setup);
23
24
25 %% validate and save setup
26 if isSetupValid(setup)
27     saveSetup('default', setup, 'override', true);
28 end
```

---

### 3.1.1 Setup file

In MATLAB a `setup` file is represented as struct. It holds the following fields

`dim`, `roi`, `coils`, and `sensors`.

Determined by the dimension of the system, `setup.dim` is either 2 or 3. The region of interest `setup.roi` is a struct with fields `x`, `y` and `z`, that holds boundary information in an array with unit meter [m]. For instance `setup.roi.x = [0,0.1]` translates to a region of interest of 10cm length in `x` direction. Yet again `setup.coils` and `setup.sensors` are arrays of structs with fields:

`coils(i).Position`, `sensors(i).Position`,  
`coils(i).Normal`, `sensors(i).Normal`,  
`coils(i).Segments`,

where  $i$  is the  $i$ -th coil/sensor. As mentioned before, spatial information (in this case the fields `Position` and `Normal`) are given as a three element vector with base unit meter [m].

The function `createEntityArray.m` is used to create new coil and sensor entries. In the current version the function takes two outer points, a normal and the number of intermediate points. Then `createEntityArray.m` returns a list of equidistant points, that represent the entity position, as well as the respective normal.

In the 3D case the coil wire is modeled by a list of points, which represent the coil shape. The function `createCoilLoop.m` creates a coil template, which has to be copied in every coil position. The function `parseCoils.m` receives a coil struct (with fields `Position` and `Normal`) and adds the `Segments` field in respect to the given coil template. Internally this uses the function `relocateStructure.m` that reorientates a given list of points.

In case a predefined voxel grid is given, the corresponding region of interest can be determined using `getROI.m`. During the creation of a new `setup` the usage of `visualizeSetup.m` is recommended.

The script in Listing 1 can be used to create the `setup` file for the 2D system shown in Figure 3. Please also examine the examples given in section 4.2, that make use of the described functionalities.

---

**Listing 2:** Create a new 2D config

---

```
1 %% create a new 2D config
2 config = [];
3
4 config.activeCoils = 1:20;
5 config.activeSensors = 1:40;
6
7 currentPattern = createCurrentPattern(config.activeCoils, 'sequential', 'current', 0.8);
8 config.currentPattern = currentPattern;
9
10 config.res = [10, 10, 1];
11
12
13 %% check compatibility with setup and save config
14 if checkCompatibility(setup, config)
15     saveConfig('configs/10.10.1/singleSequential/default', config);
16 end
```

---

### 3.1.2 Config file

config files are, similar to `setup` files, structs in MATLAB. They hold the fields

`res`, `currentPattern`, `activeCoils`, and `activeSensors`.

The simulated voxel resolution of the region of interest is given by a three element array and stored in `res`.

A setup may contain a lot of coils/sensors that are not used in a specific configuration. To optimize the simulation speed `activeCoils` and `activeSensors` are lists of coil and sensor ids, that define the subset of active setup parts.

The coil current sequence is stored in a matrix with one column for each active coil. Each row represents one pattern, where the entries refer to the applied current at the respective coil. The function `createCurrentPattern.m` provides standard patterns and multiple other presets.

The script in Listing 2 can be used to create a `config` file for the 2D system shown in Figure 3. Again please examine the examples given in Section 4.2, that make use of the described functionalities.


### 3.1.3 Setup and Config Validation

This toolbox also contain functions to check the created `setup` and `config`. The functions `isConfigValid.m` and `isSetupValid.m` validate if all required fields are set and further include checks for data inconsistencies. Additionally the function `checkCompatibility.m` checks if a given `config` is applicable for a given `setup`.

### 3.1.4 Save and Load

The `setup` and `config` datasets can be stored using `saveSetup.m` and `saveConfig.m`. These are stored using a custom `.mrxsetup` and `.mrxcfg` extension, that is based on the MATLAB internal `.mat` file format. Corresponding load functions (`loadSetup.m` and `loadConfig.m`) are available.

To store the created `setup` and `config` datasets, we recommend using the folder structure proposed in section 4.1.

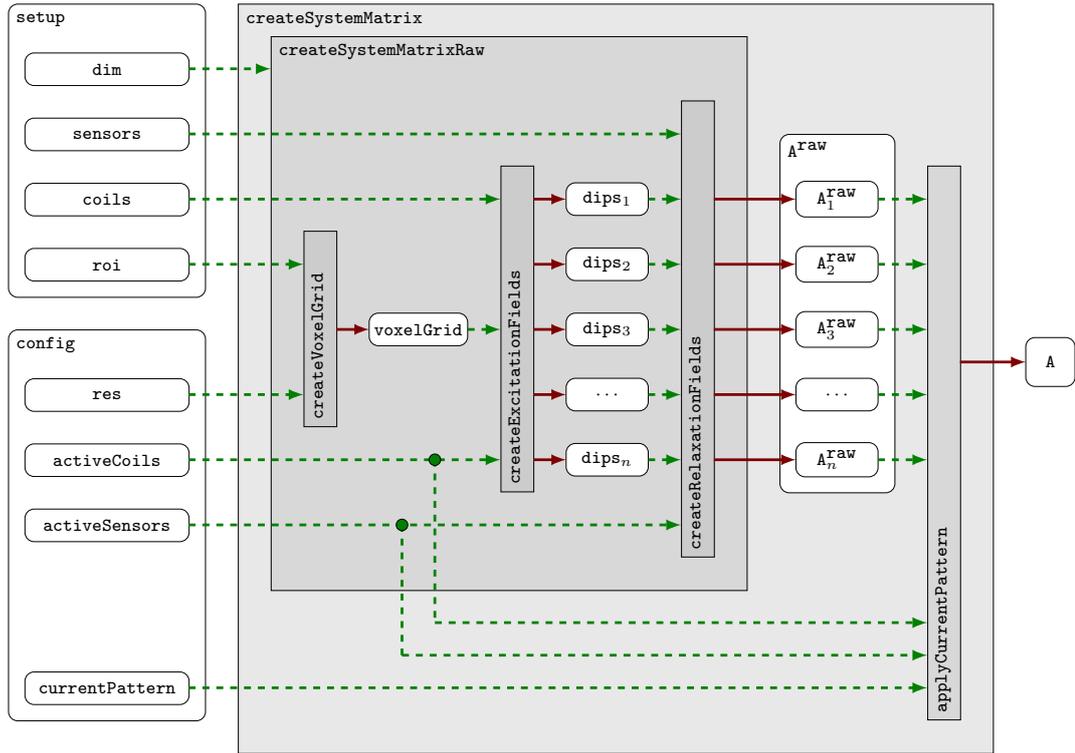


**Figure 1:** Schematics of the system matrix creation workflow. Variables are marked with a rounded corner box, functions are marked with a gray box. The variable usage is color coded: green (dashed) as input and red (solid) as output argument.

## 3.2 Simulation

The simulation module is the core element of the toolbox. All necessary files can be found in the folder `./simulation`. The main script is provided by the file `createSystemMatrix.m`. It processes a given MRX setup and configuration into a linear operator, namely a matrix  $A$ . A valid `setup` file is required, that provides information about dimension, region of interest, coil position/orientation and sensor position/orientation. Furthermore voxel/pixel resolution, lists of active coils and sensors as well as coil current patterns are provided by a `config` file. See section 3.1.3 for tools to validate the setup and config pair. Then `createSystemMatrix.m` internally calls the functions `createSystemMatrixRaw.m` and subsequently `applyCurrentPattern.m`. The internal processes, variables and function arguments are illustrated in Figure 1.

At first, `createSystemMatrixRaw.m` calls `createVoxelGrid.m`, which creates a valid `voxelGrid`. This voxel grid combines the desired resolution (`config.res`) within the demanded region of interest (`setup.roi`). Then `createExcitationFields.m` is called to calculate the magnetic fields, that are induced by an active coil with a uniform current. Internally the field calculation, hence the choice of the used mathematical model (see section 1.2), depends on the existence of `coils.Segments`. The magnetic fields are derived on the provided `voxelGrid` and then stored for each coil individually. Finally `createRelaxationFields.m` translates the fields present on the voxel grid into the magnetic response, which is acquired by the sensors (`setup.sensors`), for each coil separately in a temporary system matrix  $A_{raw}$ . The set of operators  $A_{raw}$  is finally returned by `createSystemMatrixRaw.m`.

To this point only uniform coil activation currents are used. Due to the linearity of the forward operator, we can apply any coil current sequence by using `applyCurrentPattern.m`. This derives a linear combination of the temporary system matrices  $A_{raw}$  and creates a linear operator  $A$ .

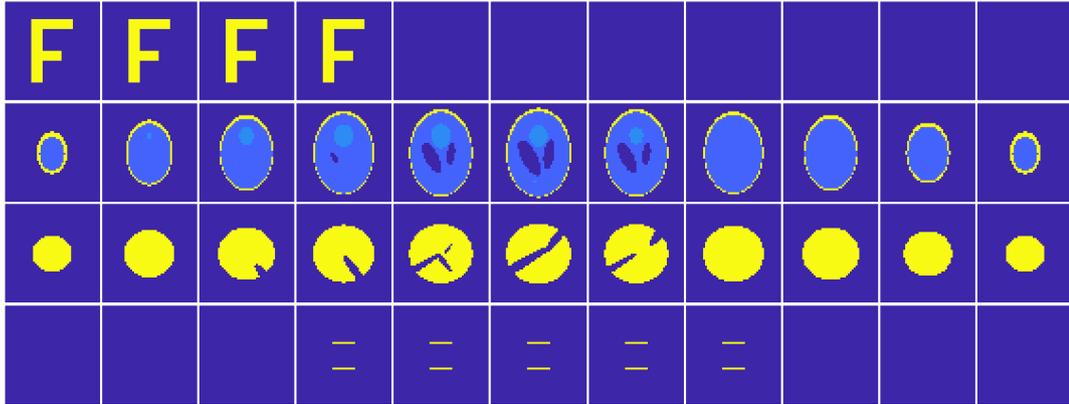

**Figure 2:** Selection of available 3D phantoms with a resolution of  $[50, 50, 15]$ . One phantom per row, each column represents a layer. Layers 1, 2, 14 and 15 are cut from illustration, due to clarity. Phantom names from top to bottom: 'F\_2', 'shepplogan3d', 'tumor' and 'fwhmdots\_0.25'.

### 3.2.1 Raw export/import

In general the simulation of the setup is the most time demanding step. It is possible to dump the dataset `ARaw` into file, which allows for flexible experimentation with varying coil currents. However this is not very suitable in case the set of active coils or sensors changes, since `ARaw` contains a fixed set of coils.

The function `exportRawSetup.m` allows for in batch export of single coil and sensor data. This allows for flexible import of any coil/sensor combination without reevaluating the system, using the function `importRawSetup.m`. Internally the simulation steps done in `exportRawSetup.m` are equivalent to those steps done by `createSystemMatrixRaw.m`. A working example (`ExampleB.m`) is presented in section 5.

## 3.3 Phantom

The function `createPhantom.m` provides multiple phantom options, including phantoms suited for reconstruction or resolution tests. `createPhantom.m` uses the function `phantom3.m`, which is a 3D reimplement of the MATLAB given phantom function and able to generate 3D phantoms that are composed of multiple ellipses. For all available phantoms, please check the help section of `createPhantom.m` and `phantom3.m`. See Figure 2 for a selection of the available phantoms.


## 3.4 Visualization

*SiMRX* provides tools to visualize the region of interest (`drawROI.m`), coils (`drawCoils.m`), sensors (`drawSensors.m`), magnetic fields (`drawFields.m`) and 3D phantoms as cuttable volume (`drawVolume.m`).

The toolbox also contains visualization functions that combine these base functions to useful tools:

The function `visualizeSetup.m` draws coils, sensors and the region of interest by passing a `setup` dataset. We recommend the usage during the creation of a new MRX setup, since it uses only the `setup` file and does not require the simulation of a given setup.

The function `compareSetups.m` is able to visually compare two setups side by side.

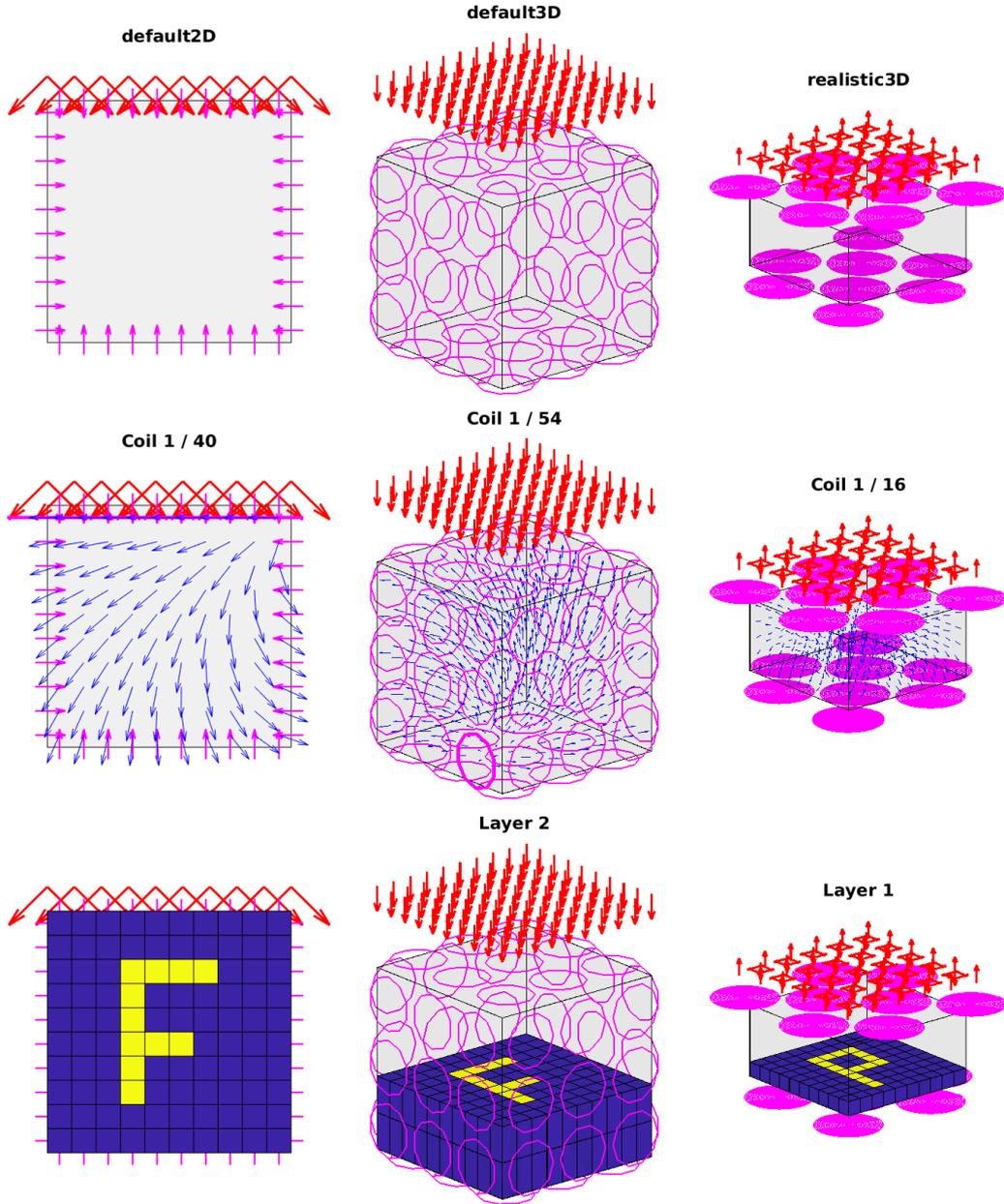

**Figure 3:** Visualization of the setups *'default2D'* (left column), *'default3D'* (center column) and *'realistic3D'* (right column). From top to bottom: visualization using `visualizeSetup.m`, `visualizeCoilFields.m` and `visualizeVolume.m`. Setups *'default2D'* and *'default3D'* are shown with phantom *'F\_2'*, setup *'realistic3D'* is shown with phantom *'P.1'*. Activation coils are colored magenta, sensors are colored red and (normalized) magnetic field lines are colored blue.

---

**Listing 3:** Proposed setup folder structure, shown for included setup '*default2D*'

---

```
setups/                                // setups base folder
├── default2D/                          // a specific setup
│   ├── README.m                       // additional information and script
│   ├── default.mrxsetup                // mrxsetup file of this setup
│   ├── configs/                       // contains all configs for this setup
│   │   ├── <resolution>/              // unique folder for each resolution
│   │   │   ├── 10.10.1/              // example resolution [10,10,1]
│   │   │   :
│   │   │   ├── <config name>/         // unique folder for each config
│   │   │   │   ├── singleSequential/ // example config 'singleSequential'
│   │   │   │   :
│   │   │   │   ├── default.mrxcfg    // mrxcfg file of this config
│   │   │   │   └── results/          // results for this setup/config/res
│   │   └── raw/                      // raw data (via exportRawSetup.m)
│   └── scripts/                      // scripts to create this config/setup
│       ├── createSetup.m              // create *.mrxsetup file
│       ├── createConfigs.m            // create *.mrxcfg in folder 'configs'
│       :
├── <my setup>/                        // my new setup
├── :
└── :
```

---

The function `visualizeCoilFields.m` provides a tool to visualize the magnetic fields induced by the activation coils. It uses the magnetic fields as an input argument, that are created by `createExcitationFields.m` during the simulation process. A scroll functionality is available to seamlessly scroll through each coil activation and current pattern respectively. A toggle between single 'coil mode' and 'current pattern mode' as well as a toggle between plotting normalized and unnormalized fields is available.

Finally `visualizeVolume.m` provides a visualization of the phantom within the corresponding setup. In the 3D case, a scroll functionality is available to cut through the volume to provide a clear view of the inner layers.

Figure 3 shows `visualizeSetup.m` and `visualizeCoilFields.m` in a 2D and 3D example. The code of the 2D example is shown in Listing 1.

## 4 Setups

This toolbox contains three setups, namely '*default2D*', '*default3D*' and '*realistic3D*'. These can be used to test the functionalities of the *SiMRX* toolbox. We will first outline the used folder structure in section 4.1 and then give a short introduction to the included setups in section 4.2.

**Please note:** Due to internal use of the provided setups for short example codes, we firmly request to not edit the setups within the *SiMRX* toolbox folder. Please create a new '*setups*' folder with copies of the provided setups outside the toolbox and make changes in there. This keeps the *SiMRX* directories clean and preserves its integrity for our internal example scripts.

### 4.1 Folder structure

The following folder structure is proposed for `setup` and its respective `config` datasets (see Listing 3 as a reference and visual representation of the folder structure).

	x [m]	y [m]	z [m]	$n_x$	$n_y$	$n_z$	SensorID	ChannelID	GroupID
Sensor 1							1		
⋮							⋮		
Sensor $N_s$							$N_s$		

(a) File: `sensors.dat`: sensor data structure

	x [m]	y [m]	z [m]
Coil 1			
⋮			
Coil $N_c$			

(b) File: `coilGrid.dat`:  
coil position data structure

	x [m]	y [m]	z [m]
Segment 1			
⋮			
Segment $N_l$			

(c) File: `coilTemplate.dat`:  
template coil data structure

	x [m]	y [m]	z [m]
Voxel 1			
⋮			
Voxel $N_v$			

(d) File: `voxelGrid.dat`:  
voxel grid data structure

	current [A]
Coil 1	
⋮	
Coil $N_c$	

(e) File: `dataset.01.currents.dat`:  
coil currents data structure

	$\Delta B$ [fT]	SensorID	ChannelID	GroupID	CoilNo
Sensor 1					1
⋮					⋮
Sensor $N_s$					1
Sensor 1					$N_c$
⋮					⋮
Sensor $N_s$					$N_c$

(f) File: `dataset.01.relax.dat`: sensor data table structure

**Table 1:** Overview of the data structures, that are used for text based datasets. **Note:** physical units are given in the brackets. The total number of sensors, coils, segments and voxel is given by  $N_s, N_c, N_l$  and  $N_v$  respectively.

For each unique setup a base folder `./setups/<my setup>` is created and all related `configs`, `scripts` and `datasets` are included in this folder. At this folders base level a `.mrxsetup` file is saved that fully describes the setup.

The folder `./setups/<setup>/configs` contains subfolders for resolution, which again contain folder for all available configs.

In `./setups/<setup>/configs/<resolution>/<config name>/` the `.mrxconfig` file is found as well as a result folder, dedicated for reconstruction results of this setup/config/resolution combination.

The folder `./setups/<setup>/raw/` is reserved for raw data exports, that has been created by the function `exportRawSetup.m` (see section 3.2.1 for the raw export/import system).

Finally the folder `./setups/<setup>/scripts/` contains scripts for this specific setup. This includes scripts for the creation of `.mrxsetup` and compatible `.mrxconfig` datasets.

## 4.2 Example setups

The respective folder `./setups/<setup>/scripts/` contain scripts for the creation of the following datasets. Furthermore the file `README.m` in `./setups/<setup>` can be executed to run all creation scripts for that setup at once.

### 4.2.1 Fully synthetic dataset: *'default2D'* and *'default3D'*

In *SiMRX* a 2D (*'default2D'*) and a 3D (*'default3D'*) example is available. These can be created using the respective MATLAB scripts `createSetup.m` and `createConfigs.m` (the 2D scrips are shown in listing 1 and 2 as seen in section 3.1.1 and 3.1.2). A visualization of both datasets is available in figure 3.4.

### 4.2.2 A 3D dataset from formatted text files: *'realistic3D'*

With *SiMRX* it is possible to load datasets from text files, save in the `.mrxsetup` and `.mrxcfg` data structure, as well as pre-processing for simulation and reconstructions tasks.

In the subfolder `./setups/realistic3D/scripts` the script `createRawDataset.m` is used to create the following files in folder `./setups/realistic3D/rawData` :

```
sensors.dat
coilGrid.dat
coilTemplate.dat
voxelGrid.dat
dataset.01.currents.dat
dataset.01.relax.dat
```

The file structure of the created data is as follows.

**Sensor Information (Table 1a, File `sensors.dat`):** The table (see Table 1a) that is used in file `sensors.dat`) stores all necessary sensor information. Each row defines a sensor unit with properties defined by the columns as follows. Columns 1 – 3 define the  $x$ ,  $y$  and  $z$  position as a translation vector. Columns 4 – 6 define the orientation of the sensors measure direction. Columns 7 – 9 are used to identify the given coil: column 7 holds an unique sensorID, column 8 holds information about the used data channel and column 9 is used to defined sensor groups.

**Coil Information (Table 1b, File `coilGrid.dat` and Table 1c, File `coilTemplate.dat`):** The file `coilGrid.dat` contains positioning information for every activation coil in the system (see Table 1b). Each row defines a coil position with columns 1 – 3 defining the translation vector for direction  $x$ ,  $y$  and  $z$ . As described in Sections 1.2 and 3.1.1, the coil is a composition of multiple conductor segments. The coils single segments are stored as a list of points in file `coilTemplate.dat` (see Table 1c).

**Voxel Grid Information (Table 1d, File: `voxelGrid.dat`):** The file `voxelGrid.dat` contains a list of used voxels in the current setup (see Table 1d). Again columns 1 – 3 define the positions of the voxel midpoints in  $x$ ,  $y$  and  $z$ .

**Current Information (Table 1e, File: dataset.01.currents.dat):** The table 1e is used as data scheme for file `dataset.01.currents.dat` and contains a list of currents in Ampere [A] that are applied to the coils. Currently, for external data, only subsequent coil patterns are supported. This means that the number of givencurrents has to fit the number of coils as of `coilGrid.dat`.

**Measurements Information (Table 1f, File: dataset.01.relax.dat):** The data structure defined in Table 1f) is used in file `dataset.01.relax.dat`. The first column contains the change in the magnetic response  $\Delta B$  [fT] after an coil activation (compare with equation (3)). Column 2 – 5 provide information about the sensorID, channelID, groupID and the coil/current pattern that was used for that measurement.

The raw dataset then found in `./setups/realistic3D/rawData` can be loaded and parsed with the scripts provided in `./setups/realistic3D/scripts`, namely `createSetup.m` and finally `createConfigs.m`. As expected, this creates a `.mrxsetup` and `.mrxcfg` file. A visualization of this dataset is available in figure 3.4. With `loadDataset.m` the data can be loaded and cleaned up in regards of faulty or unused sensors. Finally, with `simulateMeasurement.m` a simulated measurement is created based on the given dataset. Here the magnetic susceptibility  $\chi$  and amount of the particles (in [mg]) are required to define the phantom properly.

As for the other examples, a script `README.m` is available, that can be used to execute the individual scripts in the recommended order.

## 5 Examples

The folder `./exmaples` contains scripts to illustrate the workflow that is required to handle simulated and real experimental data with the *SiMRX* toolbox. It make use of the provided setups presented in section 4.2.

1. `ExampleA.m` loads the `'default2D'` setup and a eligible config file (see section 4.2.1) and simulates the corresponding system matrix.
2. `ExampleB.m` loads the `'default3D'` setup and a eligible config file (see section 4.2.1) and stores the simulated system in the raw data format (see section 3.2.1). Finally the raw data is imported and combined to a proper system matrix.
3. `ExampleC.m` loads the `'realistic3D'` setup and a eligible config file (see section 4.2.2). Then measured data are loaded, defective sensors are removed and the measured data is compared to a simulated measurement. With this example, data from real experiment can be processed.

## 6 Remarks

This simulation toolbox does not guarantee its correctness and closeness to reality. *SiMRX* is not qualified for commercial usage.

## References

- [1] Christoph Alexiou et al. “Cancer Therapy with Drug Loaded Magnetic Nanoparticles—Magnetic Drug Targeting”. In: *Journal of Magnetism and Magnetic Materials* 323.10 (2011), pp. 1404–1407 (cit. on p. 1).

- [2] Daniel Baumgarten et al. “Magnetic Nanoparticle Imaging by Means of Minimum Norm Estimates from Remanence Measurements”. In: *Medical & Biological Engineering & Computing* 46.12 (Oct. 2008), p. 1177 (cit. on p. 1).
- [3] Janic Föcke, Daniel Baumgarten, and Martin Burger. “The Inverse Problem of Magnetorelaxometry Imaging”. In: *Inverse Problems* 34.11 (Nov. 1, 2018), p. 115008 (cit. on pp. 1, 2).
- [4] James D. Hanson and Steven P. Hirshman. “Compact Expressions for the Biot–Savart Fields of a Filamentary Segment”. In: *Physics of Plasmas* 9.10 (2002), pp. 4410–4412 (cit. on p. 2).
- [5] R Hiergeist et al. “Application of Magnetite Ferrofluids for Hyperthermia”. In: *Journal of Magnetism and Magnetic Materials* 201 (Issues 1-3 1999), pp. 420–422 (cit. on p. 1).
- [6] M. Liebl et al. “Quantitative Imaging of Magnetic Nanoparticles by Magnetorelaxometry with Multiple Excitation Coils”. In: *Physics in Medicine & Biology* 59.21 (2014), p. 6607 (cit. on p. 1).