

Robust Pixel Gain Calibration with Limited Statistics

G. Blaj*, *Member, IEEE*, G. Haller, *Member, IEEE*, C. Kenney, *Member, IEEE*

I. INTRODUCTION

Pixel detectors typically display pixel-to-pixel gain variation of a few percent which result in reduced spectroscopic performance [1]. Integrating pixel detectors provide spectroscopic information directly, while photon counting detectors through threshold scans. Gain maps are used to rectify pixel-to-pixel gain variation.

Ideally, the pixel gain is calibrated with monochromatic radiation (e.g., from a synchrotron or free electron laser beamline, or laboratory x-ray sources with monochromators). However, beam time is usually valuable, and monochromatic sources are not always available [2], while x-ray tubes and radioactive sources produce complex spectra.

For calibration of ePix100a cameras, we used a Mo x-ray tube with Zr and Al filters (to optimize the Mo $K\alpha$ line) and acquired limited statistics (50 to 250 photons), sampled from a complex spectrum (where accurate fitting of a peak with charge sharing would require ≈ 1000 photons per pixel [3]).

We have developed a calibration method which relies on cross-correlating histograms of many pixel pairs and obtaining large sets of relative shifts. These were subsequently used to calculate absolute pixel shifts and corresponding pixel gains.

We demonstrate that this method yields stable gain calibration maps with an order of magnitude less statistics than required by typical approaches. Finally, we demonstrate the accuracy of the method by comparing with gain maps obtained with good statistics and monochromatic radiation at a synchrotron beamline.

The robust gain calibration method presented here can be used for any pixel detector with minimal effort and assumptions on calibration spectrum. It is particularly useful when the quantity and/or quality of calibration data is limited and repeated measurements could be difficult (e.g., previous experiments at synchrotron or free electron laser sources, FEL).

II. METHODS

We used an x-ray tube with Mo anode and Zr and Al filters to emphasize the Mo $K\alpha$ peak [3]. We acquired two data sets of ≈ 33000 frames each, with the x-ray tube operated at 45 kV and 50 kV, respectively. We collected histograms of “single pixel events” for each pixel (i.e., all 8 neighboring pixels within noise, $\leq 4\sigma$), yielding an average of 240 photons per pixel. Fig. 1(a) shows an example of pixel spectrum with

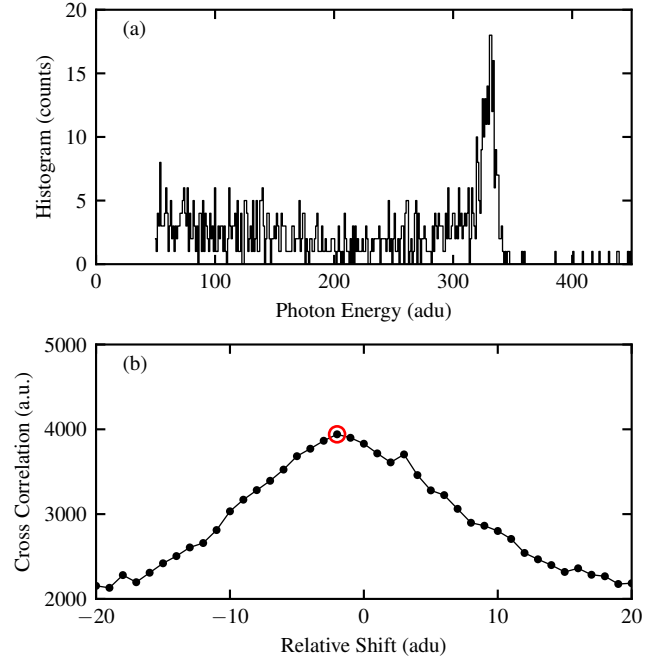


Fig. 1. Determining peak positions in statistics limited spectra is nontrivial; (a) shows a spectrum collected by one pixel; the spectrum displays significant charge sharing, a low energy shoulder, and limited statistics (on average, 240 photons per pixel near the Mo $K\alpha$ peak); (b) calculating the cross correlation of two pixel spectra and finding the maximum value yields the most likely shift between the two pixel spectra (indicated by the red circle).

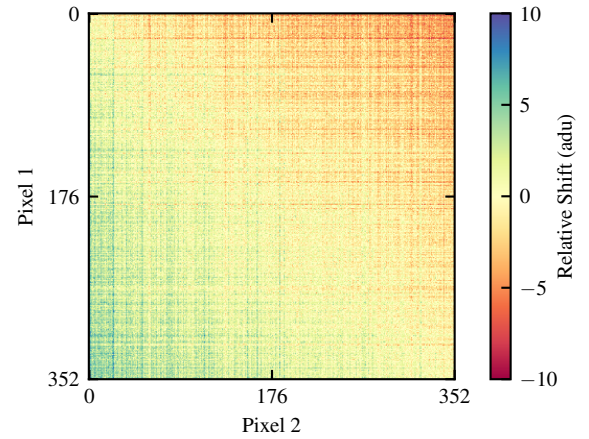


Fig. 2. The shift matrix contains relative delays for all possible pairs of pixels; individual delays could be noisy, however, averaging large numbers yields good shift estimates; we show the shift matrix for pixels in one column.

G. Blaj, G. Haller, and C. Kenney are with SLAC National Accelerator Laboratory, Menlo Park, CA 94025.

* Corresponding author: blaj@slac.stanford.edu

limited statistics. The peak position obtained through fitting

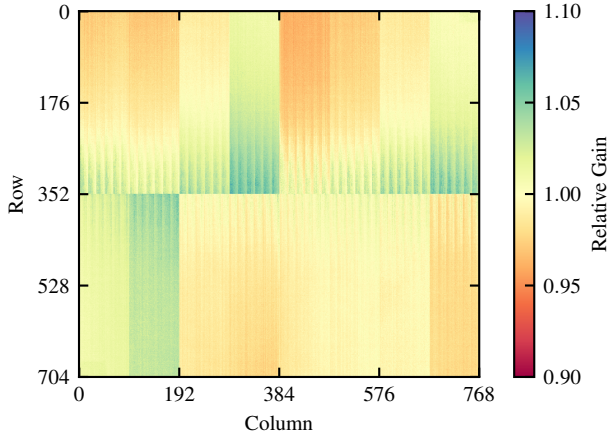


Fig. 3. Converting pixel shifts to gains yields a robust detector gain map, despite limited statistics and spectrum quality; in this ePix100a camera, 99.2% of the pixels have a relative gain in the 0.95 to 1.05 range.

or centroiding is relatively noisy [3].

A change in pixel gain results in an apparent shift of (calibration) peak position. While determining the peak shift is nontrivial with limited statistics, we developed a robust approach for spectra with a sharp transition (i.e., from a quasi-monochromatic source and/or an absorption edge filter). For each pair of pixel spectra we estimate the relative peak shift by calculating the cross correlation, finding its maximum value and corresponding delay as in Fig. 1(b). This delay is the relative shift of the two pixel spectra.

The relative shift pairs are saved in a “shift matrix” of size $N \times N$ (with N the number of pixels), see Fig. 2. With limited statistics, relative shifts of individual pairs could be imprecise, however, the column average of the shift matrix yields a robust and accurate estimate of the absolute pixel shift¹.

The algorithm complexity is $O(N^2)$; for large numbers of pixels we can optimize the run time by dividing the problem and calculating (1) shifts for pixels within each column, and (2) shifts of individual columns; their sum yields the absolute pixel i shift δ_i and the pixel gain as $g_i = 1 + \delta_i/E_0$ where E_0 is the peak position.

III. RESULTS

Fig. 3 shows the robust gain map obtained for an ePix100a camera, showing good uniformity (99.2% of the pixel gains in the 0.95 to 1.05 range). Fig. 4 shows the aggregated spectrum of all pixels, either raw (thin red line) or after gain calibration (thick black line). The spectroscopic performance is significantly improved by using a gain map [4], [5].

To test the robustness of this method, we applied it to two independent data sets obtained in different conditions. The results are shown in Fig. 5, with each black point indicating two independent gain estimations of one pixel. A linear regression yields a slope $a = 1.000 \pm 3 \times 10^{-6}$, with $R^2 = 0.987$ and residual gain noise 0.16%. This corresponds to an error due to gain calibration of 0.51 adu at the Mo K α peak, demonstrating good repeatability of the robust gain calibration with different

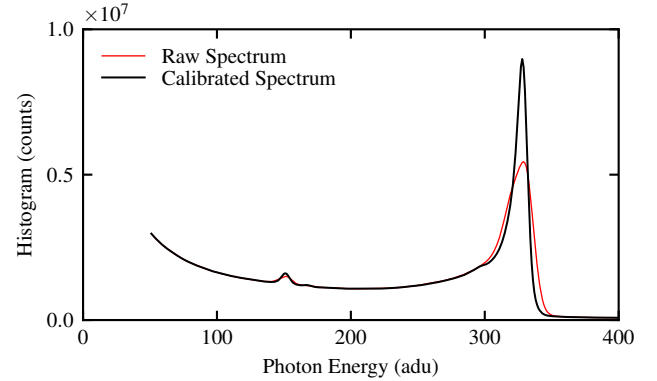


Fig. 4. Spectrum integrated over all pixels; the thin red line shows the raw spectrum and the thick black line shows the spectrum after gain calibration; the spectroscopic performance improves significantly by using a gain map.

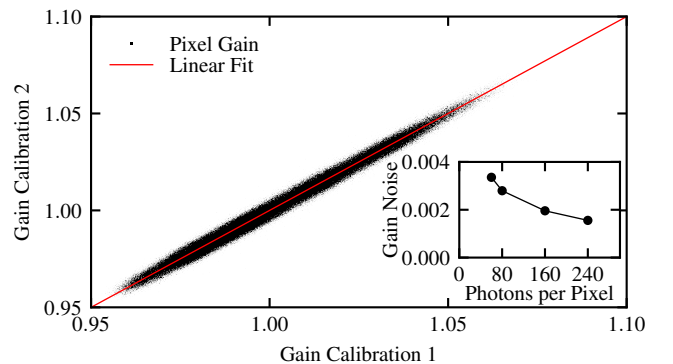


Fig. 5. Comparing gain maps obtained from two independent, statistics limited measurements; each black dot corresponds to one pixel with the two coordinates representing the two independent gains; the red line is a linear fit with a slope of $1.000 \pm 3 \times 10^{-6}$, a $R^2 = 0.987$ and residual gain noise 0.16% for the complete data set (240 photons). The inset shows the gain noise dependence on statistics; even for a reduced set with 60 photons per pixel, the gain noise is low (0.34%).

statistics limited data sets. Finally, reducing the calibration data from 240 to 60 photons per pixel results in a relatively small increase of gain noise to 0.34%.

IV. CONCLUSION

We presented a robust pixel gain calibration, demonstrating good performance and reproducibility with limited spectrum quality and statistics (an order of magnitude less than required by other methods). This method is relatively fast, with a Python program calibrating 0.5 megapixel ePix100a cameras in about 16 min on a laptop with 2.9 GHz CPU. It is also simple to implement, requiring minimal assumptions on spectrum shape.

The robust gain calibration method could be particularly useful when the quantity and/or quality of calibration data is limited and repeated measurements could be difficult (e.g., previous experimental data at synchrotron or FEL sources).

ACKNOWLEDGEMENT

Use of the Linac Coherent Light Source (LCLS) and Stanford Synchrotron Radiation Lightsource (SSRL), SLAC

¹Demonstration omitted for brevity.

National Accelerator Laboratory, is supported by the U.S. Department of Energy, Office of Science, Office of Basic Energy Sciences under Contract No. DE-AC02-76SF00515. Publication number SLAC-PUB-17459.

REFERENCES

- [1] G. Blaj, "Dead-time correction for spectroscopic photon counting pixel detectors," *Journal of Synchrotron Radiation*, 2019, in press. [Online]. Available: <https://doi.org/10.1107/S1600577519007409>
- [2] I. Kláčková, G. Blaj, P. Denes, A. Dragone, S. Göde, S. Hauf, F. Januschek, J. Joseph, and M. Kuster, "Characterization of the ePix100a and the fastccd semiconductor detectors for the European XFEL," *Journal of Instrumentation*, vol. 14, no. 01, p. C01008, 2019. [Online]. Available: <https://dx.doi.org/10.1088/1748-0221/14/01/C01008>
- [3] G. Blaj, P. Caragiulo, A. Dragone, G. Haller, J. Hasi, C. J. Kenney, M. Kwiatkowski, B. Markovic, J. Segal, and A. Tomada, "X-ray imaging with ePix100a, a high-speed, high-resolution, low-noise camera," *SPIE Proceedings*, vol. 9968, pp. 99 680J–99 680J–10, June 2016. [Online]. Available: <https://dx.doi.org/10.1117/12.2238136>
- [4] G. Blaj, A. Dragone, C. Kenney, F. Abu-Nimeh, P. Caragiulo, D. Doering, M. Kwiatkowski, B. Markovic, J. Pines, M. Weaver, S. Boutet, G. Carini, C.-E. Chang, P. Hart, J. Hasi, M. Hayes, R. Herbst, J. Koglin, K. Nakahara, J. Segal, and G. Haller, "Performance of ePix10K, a high dynamic range, gain auto-ranging pixel detector for FELs," *AIP Conference Proceedings*, vol. 2054, p. 060062, 2019. [Online]. Available: <https://doi.org/10.1063/1.5084693>
- [5] G. Blaj, D. Bhogadi, C.-E. Chang, D. Doering, C. Kenney, T. Kroll, J. Segal, D. Sokaras, and G. Haller, "Hammerhead, an ultrahigh resolution ePix camera for wavelength-dispersive spectrometers," *AIP Conference Proceedings*, vol. 2054, p. 060037, 2019. [Online]. Available: <https://doi.org/10.1063/1.5084668>