

Synchronous, Crosstalk-free Correlative AFM and Confocal Microscopies/Spectroscopies

Thales F. D. Fernandes¹, Oscar Saavedra¹, Emmanuel Margeat¹, Pierre-Emmanuel Milhiet^{1*}, and Luca Costa^{1**}

¹Centre de Biochimie Structurale (CBS), CNRS, INSERM, Univ Montpellier, 34090, Montpellier, France

*pem@cbs.cnrs.fr

**costa@cbs.cnrs.fr

ABSTRACT

Seeing is believing and consequently, microscopies have become pillars of our characterization tools to observe biological systems and assemblies. Correlative and synchronous use of different microscopies relies on the fundamental assumption of non-interference during images acquisitions. In this work, by exploring the correlative use of Atomic Force Microscopy and confocal-Fluorescence-Lifetime Imaging Microscopy (AFM-FLIM), we quantify cross-talk effects occurring during synchronous acquisition. We characterize and minimize optomechanical forces on different AFM cantilevers interfering with normal AFM operation as well as spurious luminescence from the tip and cantilever affecting time-resolved fluorescence detection. By defining non-interfering experimental imaging parameters, we show accurate real-time acquisition and two-dimensional mapping of interaction force, fluorescence lifetime and intensity characterizing morphology (AFM) and local viscosity (FLIM) of gel and fluid phases separation of supported lipid model membranes. Finally, as proof of principle by means of synchronous force and fluorescence spectroscopies, we precisely tune the lifetime of a fluorescent nanodiamond positioned on the AFM tip by controlling its distance from a metallic surface. This opens up a novel pathway of quench sensing to image soft biological samples such as membranes since it does not require tip-sample mechanical contact in contrast with conventional AFM in liquid.

Introduction

Real space observation of biomolecules at high resolution, as mentioned by Richard Feynman in his lecture at Caltech¹, is crucial to address fundamental biological questions. In this frame, several microscopy techniques have been developed during the last decades and are currently well established and available allowing characterizations at different lengths and timescales. Cryogenic Electron Microscopy (cryo-EM) offers nowadays the best spatial resolution and performs best at the molecular and sub-molecular length scales but cannot image living cells in physiological conditions². For that, a large variety of optical-based microscopies are available, historically less resolved than cryo-EM because of the diffraction limit. However, some recently developed super-resolution techniques such as Stimulated Emission Depletion Microscopy (STED)³, Stochastic Optical Reconstruction Microscopy (STORM)⁴, Photoactivated Localization Microscopy (PALM)⁵ and MINFLUX⁶ have bypassed the diffraction limit and improved the lateral resolution down to a few nanometers. While they do have some drawbacks both in terms of photo-toxicity or acquisition speed, their popularity reflects nowadays their coming of age. When it is to observe bio-interfaces at high resolution in physiological conditions, in a label-free configuration, Atomic Force Microscopy (AFM)⁷ and its fast video-rate version High-Speed AFM (HS-AFM)^{8,9} are outstanding techniques to obtain morphological and mechanical information at nano-scale. While all microscopies provide complementary physical and chemical information at different spatial and temporal resolutions, the challenge is now to get correlative observations of biological processes with two or more techniques that require important technological efforts. As examples, cryo-EM and 3D super-resolution fluorescence have been coupled to help protein localization in EM images¹⁰, as well as both STED¹¹ and STORM^{12,13} have been successfully coupled with AFM. Combining conventional fluorescence and AFM has become increasingly popular since both techniques can work in biologically relevant conditions and they acquire complementary information: specificity mapping of proteins given by specific fluorescent probes nicely complements the high-resolution AFM images. Confocal Laser Scanning Microscopy (CLSM) and AFM have been coupled at first in 1995¹⁴ to image Langmuir-Blodgett films of 10,12-pentacosadiynoic acid. In the following decades, a large number of setups following similar operational schemes based on fixed optical axes, fixed tip position, and a scanning sample were developed^{15–20}. In most setups, AFM and fluorescence images are acquired in a sequential manner, and correlative images are generated afterwards^{21,22}. However, synchronicity is of high importance as it permits to follow dynamic biological processes. In addition, it allows the acquisition of optical images while a force is applied by the AFM tip in pump and probe experiments, for instance monitoring the response of a living cell

to an AFM tip indentation or during single-molecule manipulation^{23,24}. This simultaneous modality is often not allowed due to interference and cross-talk during simultaneous image acquisition resulting in worse performance and data quality, both effects often ignored in the existing literature and only very rarely discussed^{20,25}. When coupling fluorescence (in both wide-field and confocal excitation schemes) with AFM, it is therefore of critical importance to eliminate any cross-talking between the two techniques by understanding, quantifying, and reducing the artifacts occurring during synchronous operation. We can divide artifacts into two categories:

1. AFM acquisition being affected by the excitation laser used to excite fluorescent probes. The interaction of the laser confocal spot and the AFM tip/cantilever could give rise to optomechanical/photothermal forces that can deflect the cantilever in a non-linear way which depends upon the distance from the sample. Since in time-resolved fluorescence experiments the excitation laser is continuous or pulsed at operating frequencies in the MHz regime, typically much higher than the most of cantilever resonance frequencies, the measured force (the deflection) is due to the static component of the excitation.
2. Fluorescence acquisition being affected by the presence of the AFM tip/cantilever. AFM tip luminescence gives rise to a spurious signal in the photons detection that might not be fully eliminated if there is an overlap between tip/cantilever fluorescence and the sample emission spectra.

In the case of wide-field microscopy, Total Internal Reflection excitation (TIRF) eliminates most artifacts since it limits the interaction of light only to a few hundred nanometers above the glass coverslip, where solely the tip is present, the cantilever being typically few microns axially farther. A large number of simultaneous experiments were reported in the last two decades using this modality^{26–30}. However, this type of illumination is not suitable for thick samples such as cells, since the tip will be positioned at the apical cellular membrane axial position, out of the evanescent excitation field. In a confocal scheme, the interaction of the AFM tip and the confocal spot can, unfortunately, extend to a greater distance. In this work, we focus on the coupling between a confocal spot and an AFM tip/cantilever with the aim to investigate and quantify photothermal induced deflection resulting from the presence of metal coating at the cantilever backside³¹, radiation pressure exerted on the tip³² due to the scattering of the confocal spot and finally tip/cantilever luminescence¹⁶. We show that all effects highly depend on the choice of the AFM tip and cantilever. We conclude that tip geometry, material, and metallic coating at the cantilever backside, as well as cantilever stiffness, are the key parameters to be taken into account. We also present how these common pitfalls can be avoided, or highly minimized, so there is no spurious cross-talk in the acquired images. As proof of principle we imaged synchronously the topography, lifetime, and fluorescence intensity of gel and fluid phases in Supported Lipid Bilayer (SLB), labeled with a BODIPY fluorescent probe. Finally, for the first time, to the best of our knowledge, we implemented a simultaneous AFM Force Spectroscopy (AFM-FS) and Fluorescence Spectroscopy operational scheme. In this way, we probed the force and lifetime variation due to Metal Induced Energy Transfer (MIET) between the Nitrogen-Vacancy (NV) centers in a nanodiamond attached to the apex of an AFM tip and a substrate of thin gold, therefore, enabling a novel pathway of quench sensing and opening up novel non-invasive imaging modes in liquid that can be used to image soft and fragile biological specimens.

Results

The synchronous imaging acquisition setup based on confocal excitation is described by Fig. 1. The excitation source is a supercontinuum pulsed laser where an excitation band-pass filter enables the convenient switching between different excitation wavelengths for multi-color fluorescence acquisition. The light is focused by a high numerical aperture objective. The emitted fluorescence photons are collected by the objective, focused on a pinole, recollimated and filtered and focused on an Avalanche Photodiode (APD) connected to a Time-Correlated Single Photon Counting (TCSPC) card. The AFM is mounted on top of an inverted microscope where the sample and the tip can be displaced independently with two XYZ piezo scanners. The synchronization of the AFM and the TCSPC card is made by routing the events (photons) with a Transistor-Transistor Logic (TTL) signal generated by the AFM electronics controller. Topography, intensity, and lifetime signals from the setup are used to construct simultaneously an AFM and two optical images. If operated in Force Spectroscopy mode, (FS-AFM), force, intensity, and lifetime curves versus tip-sample distance are simultaneously acquired. The lifetime was either evaluated for each image pixel fitting a decay curve with a single exponential or, in FS-AFM mode by amplitude weighting the decay curve (see Materials and Methods section). It is worth to notice that is crucial to use a short-pass filter before the APD to filter out the light used in the AFM cantilever monitoring setup³³, usually called optical beam deflection.

As explained earlier, the interaction of the focused excitation light and the AFM tip can give rise to two different phenomena: optomechanical forces and tip/cantilever luminescence. In the former case, optomechanical forces can be either due to radiation pressure (Fig. 2a) or photothermal effect (Fig. 2b). In the latter case, cantilever and/or tip luminescence (Fig. 2c) produces an additional signal creating background counts in the fluorescence channel. Radiation pressure will always be

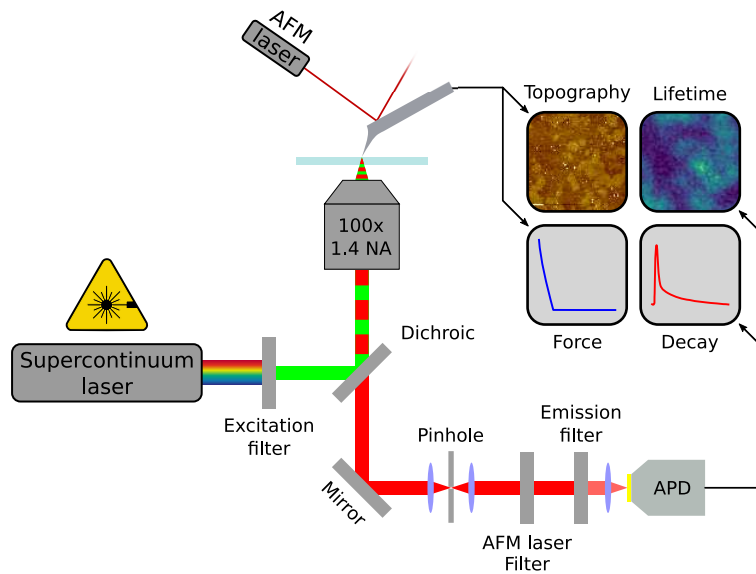


Figure 1. AFM and inverted confocal microscope setup, where the AFM and the optics are coupled by synchronizing the photons acquisition with a TTL signal from the AFM electronics controller.

present because the incoming light is scattered back (Fig. 2a) and thus a change in the momentum of light will occur. Due to momentum conservation, the tip will be deflected upwardly by a force. Photothermal deflection (Fig. 2b) is present if the cantilever is composed of two different materials, which is the case of metal-coated cantilevers, where a thin layer of metal is deposited onto the cantilever backside to improve reflectivity. Due to the different thermal expansion coefficient between the cantilever main body and the coating, when heated by the confocal spot the cantilever will be deflected. The bend will always be downwards since metals have a higher thermal expansion coefficient than silicon or silicon nitride. Therefore, the metal will expand more inducing bending to accommodate the change in length. Luminescence (Fig. 2c) can arise from two sources: 1. the material of the tip/cantilever itself can present luminescence properties, as in the case of amorphous silicon nitride, particularly important at wavelengths higher than 800 nm¹⁶. 2. The coating on the cantilever backside can be luminescent as well, as for example gold that emits in the red region of the spectrum. It is worth to notice that backscattered light will occur in all possible configurations since the tip will always scatter the incoming light of the confocal spot in all directions. However, most of the backscattered light can be easily eliminated by a notch filter when it is not sufficiently reflected by the emission filter.

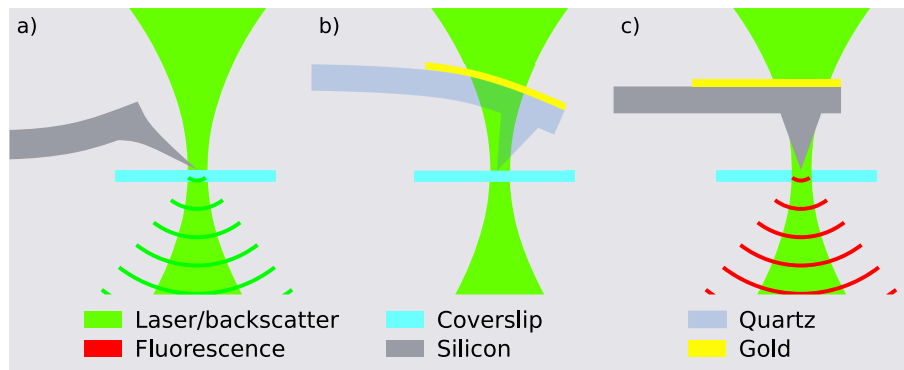


Figure 2. Possible AFM tip and confocal spot interactions: radiation pressure (a), photothermal induced deflection (b), and cantilever and tip luminescence (c).

We characterized and quantified these effects by operating the AFM in the reverse-imaging configuration³⁴. Holding the glass coverslip sample in place, the confocal spot is focused a few micrometers from the glass surface. Tip scanning is performed at constant height for different heights in open-loop. Deflection (force) and intensity, or lifetime (fluorescence), images are acquired during scan, as shown in Fig. 3. The resulting images will be of the interaction of the confocal spot with

Table 1. Optomechanical force per unit of laser power for the three AFM cantilevers

Cantilever (brand)	Optomechanical constant [pN/ μ W]
ATEC-CONT (Nanosensors)	0.740 ± 0.040
qp-BioAC (Nanosensors)	-0.672 ± 0.008
MLCT-BIO-DC (Bruker)	0.177 ± 0.005

different regions of the tip. In this operational scheme, the confocal spot is steady and operates as a “probe” while the AFM tip and cantilever are the samples to be measured. We used three AFM cantilevers with different specifications: a protruding silicon tip (Fig. 3a), a quartz-like tip with gold coating on the cantilever backside (Fig. 3b) and silicon nitride (Si_3N_4) with a gold-coated cantilever (Fig. 3c), commercially available as ATEC-CONT (Nanosensors), qp-BioAC (Nanosensors), and MLCT-BIO-DC (Bruker), respectively. The images were acquired both with the tip in focus and out of focus positions to observe the different interactions of the confocal spot (1 mW power, 532/10 nm wavelength) with the AFM tip and cantilever body. Given the large emission spectra of both silicon nitride and the metallic coating on the cantilever backside, the choice of the emission filter is not crucial if compared to fluorescence imaging of conventional organic fluorophores characterized by narrow emission spectra. In Fig. 3 we used a 690/50 nm emission filter. Panels 1, 3, and 5 show the photon intensity signal whereas panels 2, 4, and 6 display the vertical force acting on the AFM cantilever. Fig. 3a shows that the ATEC-CONT has a bright central spot in the photon intensity channel that is due to tip luminescence and presents a positive force in the range of a few hundred of piconewtons. By inserting a notch filter in the emission path, no significant change in the photon intensity channel is observed, which shows that backscattering contribution is minimal and it is completely filtered out by the band-pass emission filters. Such a positive force is a direct consequence of the radiation pressure acting on the cantilever due to the backscattered light. With the qp-BioAC (Fig. 3b) the situation is consistently different: the force is negative and it has a maximum modulus of 7 nN. The confocal spot interacts strongly with the gold reflective coating on the cantilever backside, as shown in the photon intensity images, Fig. 3b panels 1, 3, and 5. As a consequence, the interaction induces a differential heating that bends the cantilever due to photothermal effects. The photon intensity signal presents no direct correlation with the force observed and is highly dependent on the cantilever bending angle. In contrast with the previous case, the tip position is characterized by a darker spot in the photon intensity image (Fig. 3b panel 5) that is due to the fact that both excitation light and gold emission have to pass through the tip, therefore through a higher quantity of material. As a consequence, we collect less fluorescence emitted by the reflective gold coating material when the tip and confocal spot are aligned. The MLCT-BIO-DC (Fig. 3c) shows a mix of photothermal and mild radiation pressure contributions in the force signal. In this case, it is also possible to observe the interaction with the tip’s gold-coating but the thermal effects are not as drastic as before, leading to a positive total force acting on the cantilever when the tip and confocal spot are aligned.

Given the data presented in Fig. 3, the arising question is what cantilever is worth using to reduce cross-talk between AFM and time-resolved fluorescence acquisitions in correlative experiments. Fig. 4a presents time-resolved fluorescence decay curves acquired for the three AFM cantilevers of Fig. 3 with a high laser power of 1 mW to maximize signal to noise ratio for the less luminescent ones. MLCT-BIO-DC (yellow curve) has the highest photon signal by two orders of magnitude and presents a long lifetime in comparison with the qp-BioAC (red curve) and the ATEC-CONT (green curve). While ATEC-CONT and qp-BioAC will affect the measured lifetime of fluorescent probes characterized by a short (≈ 1 ns) lifetime, the MLCT-BIO-DC will affect a larger variety of fluorophores presenting even longer lifetimes. In addition, we monitored the variation of the spurious force as a function of the incident laser power (Fig. 4b). The force scales linearly with the power and the different cantilevers branch in the two categories depending on the nature of the optomechanical interaction: radiation pressure and photothermal induced deflection. It is worth to notice that the photothermal deflection computed in this case, for a laser power of 1 mW, is inferior in modulus to the 7 nN shown in Fig. 3b because in the latter case the maximum force is measured with the confocal spot focused on gold-coated cantilever backside and not on the AFM tip as in Fig. 4b. Table 1 reports the slopes from Fig. 4b presenting the total spurious force per unit of laser power.

Concluding, low laser power should be used, and the relative position between the tip and the confocal spot should be held constant, leading to a change of the cantilever equilibrium position while avoiding a change of relative deflections. In addition, our data suggest the qp-BioAC cantilevers should be employed to accurately evaluate the lifetime of the fluorescent probes within the sample (Fig. 4a). Instead, MLCT-BIO-DC cantilevers, presenting lower spurious force, are more suited for experiments where low AFM interaction forces have to be measured and in presence of a large number of fluorophores in the excitation volume, ensuring high signal to noise ratio to properly evaluate their lifetime.

The advantage of the correlative AFM-FLIM imaging mode relies on the complementarity of the two techniques. To prove that, we imaged fluid and gel phases enriched domains in biological model membranes, correlating their morphology (AFM) with their local viscosity (FLIM). We used a lipid mixture of 1,2-dipalmitoyl-sn-glycero-3-phosphocholine (DPPC), 1,2-dioleoyl-sn-glycero-3-phosphocholine (DOPC) in 1 : 1 molar ratio supplemented with 0.1% in volume of 2-(4,4-difluoro-5,7-

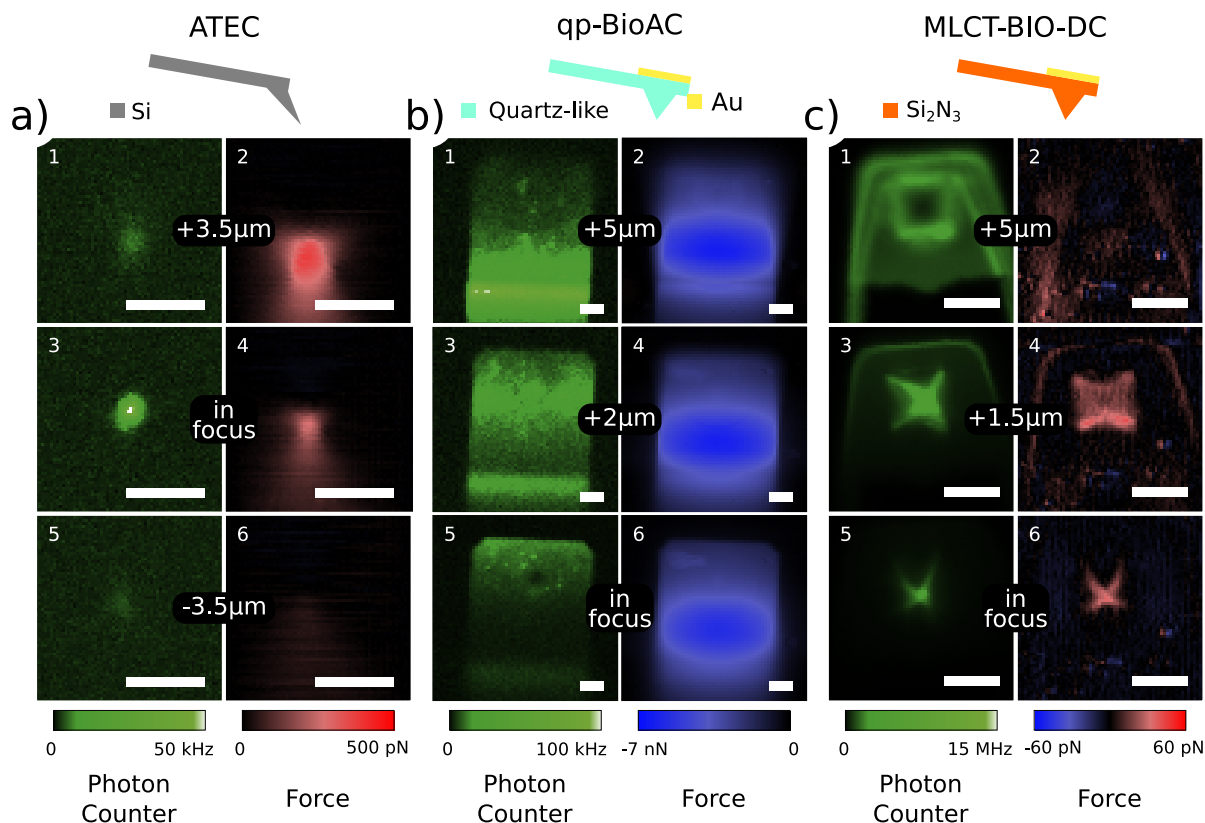


Figure 3. Simultaneous photon counts (panels 1, 3, and 5) and force (panels 2, 4, and 6) scans for three different AFM cantilevers: (a) ATEC-CONT, (b) qp-BioAC, and (c) MLCT-BIO-DC. The images were acquired using a reverse-image operational scheme where the confocal spot is a fixed probe and the AFM tip and cantilever are the measured samples. Data were collected with a 532/10 nm excitation filter, a 690/50 nm emission filter, and a laser power of 1 mW. Subfigures 1-2, 3-4, and 5-6 were acquired at different distances from the confocal spot focal plane. The scale bar is 5 μm for all images.

dimethyl-4-bora-3a,4a-diaza-s-indacene-3-dodecanoyl)-1-hexadecanoyl-sn-glycero-3-phosphocholine (BODIPY). For the simultaneous acquisition of AFM, confocal microscopy, and FLIM (shown in Fig. 5), we used the qp-BioAC cantilever and a laser power of 1 μW , ensuring, as previously explained, the highest accuracy for BODIPY lifetime evaluation while minimizing the spurious force during AFM acquisition.

Both topography (Fig. 5a) and fluorescence intensity (Fig. 5b) show that the DOPC/DPPC enriched regions (lower/higher thickness respectively) have higher/lower BODIPY concentration. Wu *et al.* showed in Giant Unilamellar Vesicles (GUVs) that higher viscosity lipids present higher lifetime, hence lifetime can be used to probe directly the lipids viscosity³⁵. Fig. 5c shows that DPPC enriched domains have higher lifetime when compared with DOPC domains, this is expected since DPPC is a gel-phase lipid, hence more viscous than DOPC, a fluid-phase lipid. The higher lifetime, above 3 ns, when compared to Wu *et al.* ≈ 1.7 ns³⁵, can be attributed to the fact that supported bilayers are sometimes affected by the substrate (glass) presence, inducing different physical behaviors in terms of lipid diffusion and potentially local viscosity when compared to free-standing bilayers in the GUVs case. The tip luminescence, dominated by low lifetime (< 1 ns) is not present in the simultaneous AFM-FLIM imaging (see section 3 of the supplementary information).

Minimizing tip/cantilever luminescence and optomechanical forces is relevant also in synchronous AFM-FS and fluorescence spectroscopy experiments. In this novel operational scheme, the AFM tip has fluorescent properties that can be correlated with the force measured by AFM while being approached or retracted from the sample. The force readout can be used to define the mechanical contact position between the sample and the fluorescent probe and, subsequently, the latter can be accurately placed with nanometric precision at a given distance from the sample and its lifetime can be probed. To demonstrate that, we show the decrease of lifetime of NV centers in a diamond nanoparticle due to MIET³⁶ while the nanodiamond is approached to a gold substrate. Fig. 6 shows an example of nanopositioning and nanomanipulation of a nanodiamond attached to the apex of an AFM tip. Fig. 6a,b show a plot of the decays for every pixel (z -position, axis), where the color scale represents the number of photons collected, acquired simultaneously with a force curve Fig. 6e,f. The experiments were

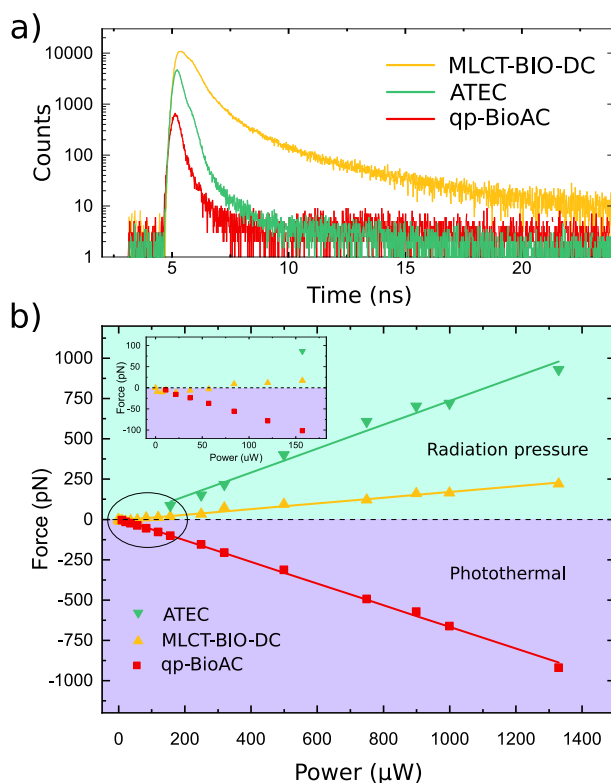


Figure 4. (a) Luminescence decay curves for different AFM cantilevers and (b) variation of the measured force as a function of the laser power. The cantilevers are qp-BioAC (red), MLCT-BIO-DC (yellow), and ATEC-CONT (green). The inset in (b) shows a zoom of the region of low power, highlighted by the ellipse. Positive force is characterized as radiation pressure while negative force is photothermal induced deflection. The slopes of the lines in (a) are: ATEC-CONT (0.74 ± 0.04 pN/μW), MLCT-BIO-DC (0.177 ± 0.005 pN/μW), and qp-BioAC (-0.672 ± 0.008 pN/μW). All curves were acquired with a 532/10 nm excitation filter, with a 690/50 nm emission filter, and 1 mW power in (a). The positive force region is colored in green to indicate the radiation pressure regime whereas the negative force region is colored in purple to indicate the photothermal regime.

performed in Dulbecco's phosphate-buffered saline (DPBS) solution. Firstly, we approached an ATEC-CONT tip holding a nanodiamond at its apex to a gold substrate, pressed with a high force (10 nN) and held in place for 2 seconds (Figs 6(a,c,e)). This allowed the nanodiamond to dislodge from the tip and to stay permanently on the gold. The tip was then retracted without the nanodiamond. This can be clearly seen from Fig 6c, where the lifetime decreases approaching the gold and never recovers back during tip retraction from gold. In a second experiment, showed in Figs 6(b,d,f), we used a qp-BioAC holding a stable nanodiamond at the tip's apex (see Materials and Methods section). In this configuration, the fluorescent nanodiamond is approached, pressed with lower force (5 nN), and retracted from the gold substrate without being dislodged from the AFM tip. The lifetime and force curves are symmetric, Figs 6(d,f), showing that the interaction is reversible. The quenching in lifetime when approaching gold, or any metal, is expected due to MIET³⁶. The radiation rate will change due to the coupling of the nanodiamond electromagnetic near-field with the plasmons of the gold. In the vertical position (x -axis) of Figs. 6(e,f) the approach and retract curves have been “unfolded”, where approach/retract have negative/positive position values. A folded version of the force curves is presented in section 2 of the supplementary information together with the lifetime curve.

Discussion

The mutual interaction of the confocal spot and the AFM tip/cantilever can give rise to a cross-talk in the AFM and/or fluorescence images. Two main forces have been explored: radiation pressure and photothermal induced deflections. Miranda *et al.* have calculated the deformation due to laser heating on gold-coated cantilevers and found displacements of 9 nm and forces around 3 nN with 6 times less power than used in Fig. 3²⁵. Radiation pressure forces can be estimated to be in the range of few tens/hundreds of piconewtons, in agreement with our results (see section 1 of the supplementary information). While in standard confocal imaging a laser power inferior to 100 μW is usually employed, in correlative STED-AFM¹² a much higher

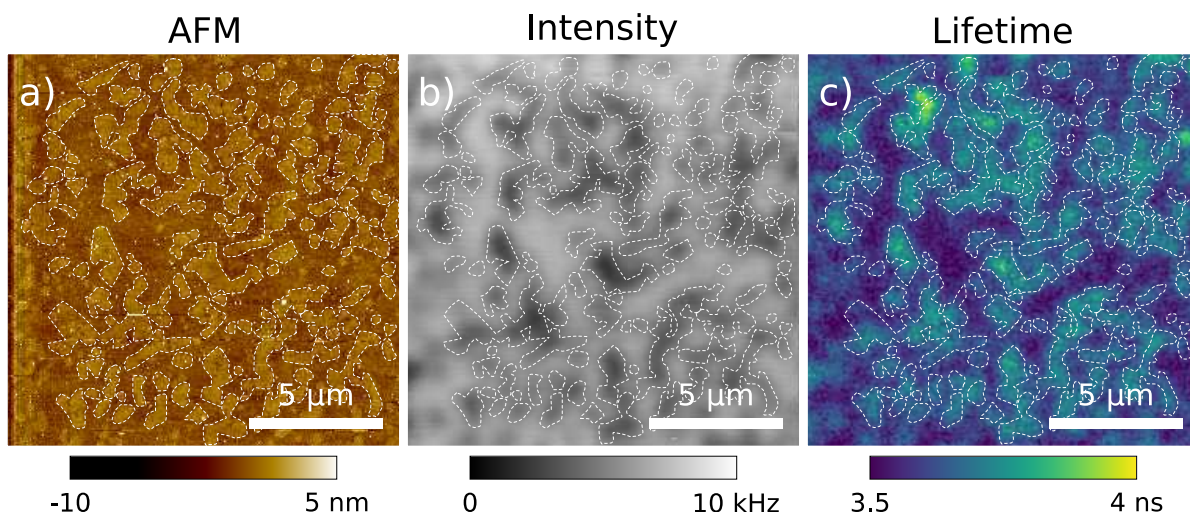


Figure 5. AFM morphology (a), fluorescence intensity (b), and lifetime (c) of DOPC:DPPC (1 : 1 molar ratio) labeled with 0.1% BODIPY model membranes. The three images were synchronously acquired. The scale bar is 5 μm , the power used was 1 μW , with a 488/10 nm excitation filter and a 525/39 nm emission filter. The tip used was a qp-BioAC. The dashed white line serves as visual aid for the separation of the two domains.

power in the mW range is required, suggesting that simultaneous STED and AFM will produce optomechanical forces from several hundreds of pN to the nN range (Fig. 4b), causing large tip deflections when using soft AFM cantilevers (< 0.1 N/m stiffnesses). The inset of Fig. 4b, however, suggests that even at reasonably low laser power (10-50 μW) optomechanical forces are present and can be quantified in few tens of pN. The same applies to the overall attractive force in the case of qp-BioAC. Those forces are not negligible and should be taken into account when measuring by AFM particularly soft interfaces such as living cells plasma membrane³⁷ or when unfolding single molecules³⁸. Moreover, the optomechanical interaction will not be constant if the alignment between the excitation laser and the AFM tip is not kept fixed during data acquisition, for instance, if the confocal spot is scanned and the tip is fixed as in the case of AFM coupled with a confocal microscope with a spinning disk operational scheme²³ or with STED. While MLCT-BIO-DC minimizes optomechanical forces, quartz-like cantilevers such as qp-BioAC perform better for time-resolved fluorescence measurements due to lower tip luminescence. We have shown that using a low excitation power of 1 μW (more generally inferior to 10 μW) it is possible to image synchronously model membrane morphology, fluorescence intensity and lifetime minimizing optomechanical spurious forces (< 10 pN) and tip and gold coating on cantilever backside luminescence, therefore minimizing any tip-related artifact.

Concerning the nanopositioning experiment presented in Fig. 6, the lifetime decrease occurs in absence of mechanical contact between the fluorescent probe and the gold substrate, at short distances where the AFM is not sensitive enough to detect any interaction force. In perspective, this enables the possibility to obtain topographical images in liquid acquired in absence of AFM tip-sample mechanical contact. It represents an important advance when it is to image soft, fluorescent-labeled, samples that are largely deformed by conventional AFM due to their extreme softness as it occurs in the case of living cell membranes. In our case, only a small drop (8-5%) in nanodiamond lifetime is observed in Fig. 6c,d, respectively. Tisler *et al.* have observed a 35% MIET efficiency when approaching a single NV center 25 nm nanodiamond attached to an AFM tip on top of graphene, a material that behaves like a metal¹⁸. In our case, the small MIET efficiency observed in Fig. 6 could be attributed to either the larger nanodiamond size (≈ 40 nm) or the more abundant number of NV centers (≈ 15) that are randomly oriented in the nanodiamond core. In the case of the experiment shown by Tisler *et al.*¹⁸ a complete lifetime spectroscopical change in respect to the nanodiamond-sample distance is not reported. Indeed, in their case, the lifetime variation was measured only when the nanodiamond was in mechanical contact with graphene or the underlining substrate. We stress that insight in the near field non-contact interaction should be probed through the acquisition of the full lifetime approach curve where a fluorescence decay is recorded for each z-position of the tip/emitter, as it has been shown by Buchler *et al.* in the case of a gold nanoparticle approached to a glass coverslip, validating the theory predicted by Drexhage³⁶ in 1970. In addition, we here stress that the full force approach curve should be simultaneously acquired to the lifetime curve for each z-position (Fig. 6). In this way, the precise mechanical contact point can be identified from the AFM interaction force and consequently an accurate fluorescent tip-sample distance can be measured.

In summary, our results impose important experimental constraints for synchronous operation and interference-free correlative AFM and confocal time-resolved fluorescence. In this work, we have quantified in fraction of pN/ μW the radiation

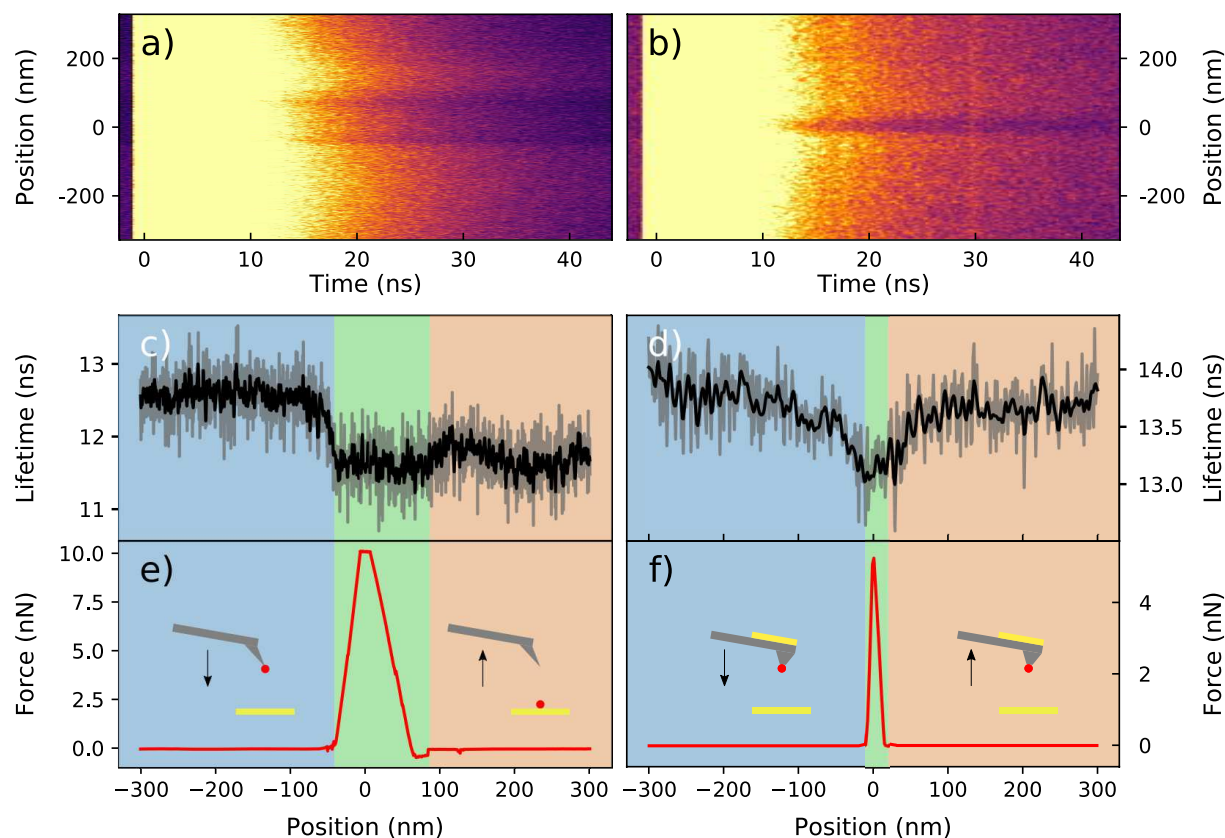


Figure 6. Simultaneous force curve and nanodiamond lifetime measurements approaching and retracting a gold substrate. Photons are synchronized with TTL from the AFM that enables the acquisition of decay curves for every z -position of the tip. (a,b) show tip z -position versus time (decay curves). (c,d) show lifetime of the nanodiamond in function of the tip position. (e,f) show standard force curves unfolded so approach and retract curves are separated. The blue region in (c-f) consists of the approach curve whereas red region the retract curve. The green region highlights the region where the tip is in contact with the sample. A 532/10 nm excitation filter was used with a 709/167 nm emission filter with power of 300 μ W. The light gray curve in (c,d) is the curve without binning (filter) whereas in the black curve a spatial binning of 9 neighbors and 5 temporal neighbors was used. Spatial binning is also applied to (a,b) to aid visualization. An ATEC-CONT was used in (a,c,f) whereas in (b,d,f) a qp-BioAC with a plateau was used. Tip schematic in (e,f) show an AFM tip with a nanodiamond (in red) and the substrate (yellow). In the retract curve of (e) the nanodiamond is left on the surface whereas in (f) it is kept at the tip's end.

pressure and photothermal detection forces acting on AFM cantilevers composed of different materials and in the case of different tip geometries. While quartz-like AFM probes, minimizing tip luminescence, are more suited when fluorescence lifetime has to be measured, the overall optomechanical forces are minimized if radiation pressure and photothermal effect contributions are comparable, canceling each other. However, our findings allow for the synchronous collection of AFM and FLIM images as well as for the simultaneous correlation of Force Spectroscopy (AFM-FS) with fluorescence lifetime. We have shown the correlation of the topography (AFM) with the local viscosity (FLIM) in the presence of gel and fluid phase separated domains in model biological membranes labeled with a BODIPY fluorophore. Finally, we have shown that the simultaneous correlation of force and fluorescence lifetime reveals a unique potential that we have used to reversibly tune the spectroscopic properties of a fluorescent nanodiamond attached to an AFM tip once it is placed in close vicinity with a gold substrate.

Materials and Methods

Atomic Force and Confocal Microscopies images were acquired using a Nanowizard 4 (JPK Instruments, Bruker) equipped with a Tip Assisted Optics (TAO) module and a Vortis-SPM control unit. The AFM cantilever detection system employs an infrared low-coherence light source, with emission centered at 980 nm and further long-pass filtered at 850 nm with a Schott RG850 filter. The AFM head was mounted on a Zeiss inverted optical microscope. A custom-made confocal microscope was coupled to the AFM using a Rock-PP supercontinuum laser (Leukos) as an excitation source operating at 20 MHz (50 ns pulses of width < 100 ps). We used an oil immersion objective with a 1.4 numerical aperture (Plan-Apochromat 100x, FWD=0.17mm, Zeiss) and a pinhole of 100 μm diameter size (P100D, Thorlabs). We used an avalanche photodetector (SPCM-AQR, PerkinElmer) connected to an SPC-150 (Becker&Hickl) TCSPC card to collect fluorescence signals. The light source of the AFM detection system was filtered out in the fluorescence emission path using an ET800sp short pass filter (Chroma). The excitation laser power was quantified with an S121C detector, located in the excitation path, and a PM100 energy meter interface (both purchased from Thorlabs). We used 532/10 nm or 488/10 nm excitation filters and 709/167 nm or 525/39 nm (from Semrock and Thorlabs, respectively) as emission filters.

During simultaneous AFM/FLIM images the AFM tip and confocal spot positions were fixed by co-aligning them. Full alignment is obtained in two steps: 1. The tip is at first roughly aligned using white field illumination and then the alignment is fine-tuned using the increase of tip luminescence. It is worth to notice that when tip and confocal spot are aligned, some AFM tips, such as qp-BioAC (Nanosensors), show colocalized tip luminescence and optical forces maxima, whereas so-called “protruding AFM tips” do not show maxima colocalization. 2. The sample is imaged and any mismatch between topography and FLIM is adjusted by correcting the tip position.

Force/Lifetime curves. Acquisition of approach-retract curves with simultaneous collection of force, photon intensity and fluorescent nanodiamond lifetime was performed using 532/10 nm as excitation filter and 709/167 nm as emission filter with a laser excitation power of 300 μW . AFM force spectroscopy was performed in static mode at a speed of 5 nm/s and a sampling rate of 20 Hz. Lifetime was evaluated by amplitude weighting the decay curve for times greater than 6 ns. All measurements were performed at room temperature and under liquid (DPBS) conditions. The tip was held constant and the sample was approached to the tip using the axial sample piezo of the TAO module. The tip-confocal alignment was performed similarly to the method described above with a key change: we focused the confocal spot when the tip is in mechanical contact with the sample surface, hence the maximum signal from the nanodiamond/gold happens at the contact point (Fig. 6).

FLIM of model membranes. Images were obtained using a 488/10 nm excitation filter and a 535/39 nm emission filter with a laser excitation power of 1 μW . AFM images of the membrane were acquired in contact mode using a rectangular qp-BioAC cantilever with a stiffness of 0.3 N/m at constant force of 1.5 nN at a line imaging rate of 0.25 Hz with 128 lines \times 128 pixels and with a scan size of 15 μm \times 15 μm . Images were collected with a sample scan with the AFM tip and excitation beam kept aligned and steady. Lifetime was fitted with SPCimage (Becker&Hickl) using a single exponential and spatial binning of 1.

AFM tips. qp-BioAC and ATEC-CONT were purchased from Nanosensors and MLCT-BIO-DC was purchased from Bruker. Custom tips, presenting a plateau at the tip apex, were purchased from Nanoworld: they were fabricated modifying a qp-BioAC cantilever. An electron beam deposited carbon tip was grown on top of a thin layer of gold. Subsequently, a flat circular plateau of 50 nm diameter was obtained at the tip apex by means of Field Ion Beam (FIB) to aid fluorescent nanodiamond fixation and stability.

Tip calibration. In all AFM experiments, the inverse optical lever sensitivity and lever stiffness of qp-BioAC and ATEC-CONT cantilevers were calibrated using a combination of a Sader³⁹ and thermal⁴⁰ methods (called “contact-free” method in the JPK AFM instruments). MLCT-BIO-DC triangular cantilever optical lever sensitivity was evaluated with the acquisition of an approach-retract curve (deflection vs z scanner displacement curve) on a glass rigid surface and cantilevers stiffness were quantified consequently with a thermal method.

Reverse-tip imaging. Simultaneous photon counts and force maps presented in Fig. 3 were obtained by scanning the tip both in focus and out of focus positions at a constant height in AFM open loop, therefore in the absence of active feedback on the cantilever deflection and with the excitation beam kept steady in space during data collection. Images were acquired using a 532/10 nm excitation filter, a measured power of 1 mW and an emission filter of 690/50 nm using 64 lines \times 64 pixels, a line rate of 6 Hz and scan size of 10 μm \times 10 μm for ATEC-CONT probes, 15 μm \times 15 μm for the MLCT-BIO-DC and 30 μm \times 30 μm for the qp-BioAC.

Model membranes. 1,2-dipalmitoyl-sn-glycero-3-phosphocholine (DPPC) and 1,2-dioleoyl-sn-glycero-3-phosphocholine (DOPC) were purchased from Avanti Polar Lipids Inc. 2-(4,4-difluoro-5,7-dimethyl-4-bora-3a,4a-diaza-s-indacene-3-dodecanoyl)-1-hexadecanoyl-sn-glycero-3-phosphocholine (BODIPY) was purchased from Thermofisher. Chloroform and methanol were purchased from Sigma-Aldrich. Membrane experiments were performed in DPBS buffer solution (GIBCO) filtered before use with an inorganic membrane filter (0.22 μm pore size, Whatman International Ltd). DPPC and DOPC were individually dissolved in chloroform:methanol (v:v 1:1) and supplemented with 0.1% BODIPY. The solvent was then evaporated to dryness under nitrogen flow to obtain a thin film spread in a glass tube. The dried lipids films were hydrated with DPBS buffer solution, previously heated at 60°C, until a final total concentration of 0.2 mM. The tube was later subjected to cycles of vortex mixing and heating at 60°C. The vesicle suspension was extruded with a polycarbonate membrane filter (100 nm pore size, Whatman, purchased from Avanti Lipids). Circular glass coverslips (2.8 cm diameter, 165 μm thick, purchased from Marienfeld) were cleaned by a cycle of sonication in KOH for 15 min, followed by a second cycle of sonication in deionized water for 15 min. Then, the glass coverslips were exposed to plasma (Expanded Plasma Cleaner PDC-002, Harrick Scientific Corporation) at high RF power level for 15 min. Supported lipid bilayers (SLBs) were obtained by vesicles fusion method^{41–43}. 300 μL of vesicles suspensions were deposited onto cleaned glass coverslips, previously mechanically fixed in an Attofluor chamber (Thermofisher), and incubated for 30 min at 70°C. Afterwards, the samples were rinsed several times with buffer solution to avoid unfused vesicles, always keeping the substrates hydrated and imaged by correlative confocal-AFM after 24 hours.

Nanodiamonds and gold substrate. Circular glass coverslips (Marienfeld) were subjected to sonication cycles in acetone (5 minutes) and ethanol (5 minutes), and subsequently dried under nitrogen flow. 5 nm or 10 nm of gold were evaporated on top of coverslips by thermal evaporation method at a rate of 1 Å/s (AS053 Oerlikon thermal evaporating source). Fluorescent Nanodiamonds (FNDs) were purchased from FND Biotech (Taiwan). A solution of 40 nm diameter nanodiamonds, each containing in average 15 NV centers, at a concentration of 0.1 mg/mL was sonicated for 15 min and subsequently diluted to 0.1 $\mu\text{g/mL}$. 5 μL were deposited on both gold-coated and glass coverslips and left to dry. Finally, nanodiamonds were imaged by correlative confocal and AFM in deionized DPBS solution. Nanodiamond grafting at the AFM tip apex was performed on both ATEC-CONT and qp-BioAC tips. The ATEC-CONT was scanned on a sample of gold (5 nm thick) with previously deposited nanodiamonds on top. A single nanodiamond was first localized through the acquisition of an AFM image and then it was picked from the surface by pressing on it with a few nanonewtons³⁴. The modified qp-BioAC, with a circular plateau at the tip apex, was instead incubated for 15 min in poly-L-lysine (PLL). Then it was used to localize a single nanodiamond previously deposited on a glass coverslip and picked it upon exertion of a force of 1 nN. The use of PLL permanently sets the nanodiamond to be at the tip apex if low forces are used to avoid dislodging the nanodiamond.

Acknowledgments

L.C. and T.F. acknowledge support from CNRS Momentum program (2017). O.S. acknowledges funding from the European Union's Horizon 2020 research and innovation program under the Marie Skłodowska-Curie grant agreement No. 721874 (SPM 2.0). The CBS is a member of the France-BioImaging (FBI) and the French Infrastructure for the Integrated Biology (FRISBI), 2 national infrastructures supported by the French National Research Agency (ANR-10-INBS-04-01 and ANR-10-INBS-05 respectively). The authors acknowledge Francesco Pedaci for discussions and Martina Di Muzio and the Micro-FabSpace and Microscopy Characterization Facility, Unit 7 of ICTS "NANBIOSIS" from CIBER at IBEC, for providing the gold-coated coverslips.

Author contributions

E.M., P.-E. M., L.C. and T.F.D.F conceived the research. T.F.D.F and L.C. conducted the experiments. T.F.D.F analyzed the data and conceived the Python code for data analysis. T.F.D.F, O.S., L.C. and E.M. set up the instrument. All authors reviewed the manuscript.

References

1. Feynman, R. P. There's plenty of room at the bottom. *California Inst. Technol. Eng. Sci. magazine* (1960).

2. Cheng, Y. Single-particle cryo-EM—How did it get here and where will it go. *Science* **361**, 876–880, DOI: [10.1126/science.aat4346](https://doi.org/10.1126/science.aat4346) (2018).
3. Willig, K. I., Rizzoli, S. O., Westphal, V., Jahn, R. & Hell, S. W. STED microscopy reveals that synaptotagmin remains clustered after synaptic vesicle exocytosis. *Nature* **440**, 935–939, DOI: [10.1038/nature04592](https://doi.org/10.1038/nature04592) (2006).
4. Rust, M. J., Bates, M. & Zhuang, X. Sub-diffraction-limit imaging by stochastic optical reconstruction microscopy (STORM). *Nat. Methods* **3**, 793–796, DOI: [10.1038/nmeth929](https://doi.org/10.1038/nmeth929) (2006).
5. Manley, S. *et al.* High-density mapping of single-molecule trajectories with photoactivated localization microscopy. *Nat. Methods* **5**, 155–157, DOI: [10.1038/nmeth.1176](https://doi.org/10.1038/nmeth.1176) (2008).
6. Balzarotti, F. *et al.* Nanometer resolution imaging and tracking of fluorescent molecules with minimal photon fluxes. *Science* **355**, 606–612 (2017).
7. Dufrène, Y. F. *et al.* Imaging modes of atomic force microscopy for application in molecular and cell biology. *Nat. Nanotechnol.* **12**, 295–307, DOI: [10.1038/nnano.2017.45](https://doi.org/10.1038/nnano.2017.45) (2017).
8. Ando, T. High-speed atomic force microscopy coming of age. *Nanotechnology* **23**, 062001, DOI: [10.1088/0957-4484/23/6/062001](https://doi.org/10.1088/0957-4484/23/6/062001) (2012).
9. Nasrallah, H. *et al.* Imaging artificial membranes using high-speed atomic force microscopy. In *Atomic Force Microscopy*, 45–59 (Springer, 2019).
10. Kopek, B. G., Shtengel, G., Xu, C. S., Clayton, D. A. & Hess, H. F. Correlative 3D superresolution fluorescence and electron microscopy reveal the relationship of mitochondrial nucleoids to membranes. *Proc. Natl. Acad. Sci.* **109**, 6136–6141, DOI: [10.1073/pnas.1121558109](https://doi.org/10.1073/pnas.1121558109) (2012).
11. Harke, B., Chacko, J. V., Haschke, H., Canale, C. & Diaspro, A. A novel nanoscopic tool by combining AFM with STED microscopy. *Opt. Nanoscopy* **1**, 3, DOI: [10.1186/2192-2853-1-3](https://doi.org/10.1186/2192-2853-1-3) (2012).
12. Odermatt, P. D. *et al.* High-Resolution Correlative Microscopy: Bridging the Gap between Single Molecule Localization Microscopy and Atomic Force Microscopy. *Nano Lett.* **15**, 4896–4904, DOI: [10.1021/acs.nanolett.5b00572](https://doi.org/10.1021/acs.nanolett.5b00572) (2015).
13. Dahmane, S. *et al.* Nanoscale organization of tetraspanins during HIV-1 budding by correlative dSTORM/AFM. *Nanoscale* **11**, 6036–6044, DOI: [10.1039/C8NR07269H](https://doi.org/10.1039/C8NR07269H) (2019).
14. Hillner, P., Walters, D., Lal, R., Hansma, H. & Hansma, R. Combined Atomic Force and Confocal Laser Scanning Microscope. *Microsc. Microanal.* **1**, 127–130, DOI: [10.1017/S1431927695111277](https://doi.org/10.1017/S1431927695111277) (1995).
15. Gradinaru, C. C., Martinsson, P., Aartsma, T. J. & Schmidt, T. Simultaneous atomic-force and two-photon fluorescence imaging of biological specimens in vivo. *Ultramicroscopy* **99**, 235–245, DOI: [10.1016/j.ultramic.2003.12.009](https://doi.org/10.1016/j.ultramic.2003.12.009) (2004).
16. Kassies, R. *et al.* Combined AFM and confocal fluorescence microscope for applications in bio-nanotechnology. *J. Microsc.* **217**, 109–116, DOI: [10.1111/j.0022-2720.2005.01428.x](https://doi.org/10.1111/j.0022-2720.2005.01428.x) (2005).
17. Bek, A. *et al.* Fluorescence Enhancement in Hot Spots of AFM-Designed Gold Nanoparticle Sandwiches. *Nano Lett.* **8**, 485–490, DOI: [10.1021/nl072602n](https://doi.org/10.1021/nl072602n) (2008).
18. Tisler, J. *et al.* Single Defect Center Scanning Near-Field Optical Microscopy on Graphene. *Nano Lett.* **13**, 3152–3156, DOI: [10.1021/nl401129m](https://doi.org/10.1021/nl401129m) (2013).
19. Tetienne, J.-p. *Un microscope de champ magnétique basé sur le défaut azote-lacune du diamant: réalisation et application à l'étude de couches ferromagnétiques ultraminces*. Ph.D. thesis, École normale supérieure de Cachan - ENS Cachan (2014).
20. He, Y., Govind Rao, V., Cao, J. & Lu, H. P. Simultaneous spectroscopic and topographic imaging of single-molecule interfacial electron-transfer reactivity and local nanoscale environment. *The journal physical chemistry letters* **7**, 2221–2227, DOI: [10.1021/acs.jpcllett.6b00862](https://doi.org/10.1021/acs.jpcllett.6b00862) (2016).
21. Meller, K. & Theiss, C. Atomic force microscopy and confocal laser scanning microscopy on the cytoskeleton of permeabilised and embedded cells. *Ultramicroscopy* **106**, 320–325, DOI: [10.1016/j.ultramic.2005.10.003](https://doi.org/10.1016/j.ultramic.2005.10.003) (2006).
22. Laskowski, P. R. *et al.* High-Resolution Imaging and Multiparametric Characterization of Native Membranes by Combining Confocal Microscopy and an Atomic Force Microscopy-Based Toolbox. *ACS Nano* **11**, 8292–8301, DOI: [10.1021/acs.nano.7b03456](https://doi.org/10.1021/acs.nano.7b03456) (2017).
23. Efremov, Y. M. *et al.* Anisotropy vs isotropy in living cell indentation with AFM. *Sci. Reports* **9**, 5757, DOI: [10.1038/s41598-019-42077-1](https://doi.org/10.1038/s41598-019-42077-1) (2019).

24. Beicker, K., O'Brien, E. T., Falvo, M. R. & Superfine, R. Vertical light sheet enhanced side-view imaging for afm cell mechanics studies. *Sci. Reports* **8**, 1504– (2018).
25. Miranda, A., Martins, M. & De Beule, P. A. A. Simultaneous differential spinning disk fluorescence optical sectioning microscopy and nanomechanical mapping atomic force microscopy. *Rev. Sci. Instruments* **86**, 093705, DOI: [10.1063/1.4931064](https://doi.org/10.1063/1.4931064) (2015).
26. Sarkar, A., Robertson, R. B. & Fernandez, J. M. Simultaneous atomic force microscope and fluorescence measurements of protein unfolding using a calibrated evanescent wave. *Proc. Natl. Acad. Sci.* **101**, 12882–12886, DOI: [10.1073/pnas.0403534101](https://doi.org/10.1073/pnas.0403534101) (2004).
27. Gump, H., Stahl, S. W., Strackharn, M., Puchner, E. M. & Gaub, H. E. Ultrastable combined atomic force and total internal fluorescence microscope. *Rev. Sci. Instruments* **80**, 063704, DOI: [10.1063/1.3148224](https://doi.org/10.1063/1.3148224) (2009).
28. Ortega-Esteban, A. *et al.* Fluorescence Tracking of Genome Release during Mechanical Unpacking of Single Viruses. *ACS Nano* **9**, 10571–10579, DOI: [10.1021/acsnano.5b03020](https://doi.org/10.1021/acsnano.5b03020) (2015).
29. Uchihashi, T., Watanabe, H., Fukuda, S., Shibata, M. & Ando, T. Functional extension of high-speed AFM for wider biological applications. *Ultramicroscopy* **160**, 182–196, DOI: [10.1016/j.ultramic.2015.10.017](https://doi.org/10.1016/j.ultramic.2015.10.017) (2016).
30. Umakoshi, T., Fukuda, S., Iino, R., Uchihashi, T. & Ando, T. High-speed near-field fluorescence microscopy combined with high-speed atomic force microscopy for biological studies. *Biochimica et Biophys. Acta (BBA) - Gen. Subj.* **1864**, 0–1, DOI: [10.1016/j.bbagen.2019.03.011](https://doi.org/10.1016/j.bbagen.2019.03.011) (2020).
31. Ramos, D., Tamayo, J., Mertens, J. & Calleja, M. Photothermal excitation of microcantilevers in liquids. *J. Appl. Phys.* **99**, 124904, DOI: [10.1063/1.2205409](https://doi.org/10.1063/1.2205409) (2006).
32. Favero, I. & Karrai, K. Optomechanics of deformable optical cavities. *Nat. Photonics* **3**, 201–205, DOI: [10.1038/nphoton.2009.42](https://doi.org/10.1038/nphoton.2009.42) (2009).
33. Araghi, H. Y. & Paige, M. F. Insight into diacetylene photopolymerization in Langmuir-Blodgett films using simultaneous AFM and fluorescence microscopy imaging. *Surf. Interface Analysis* **49**, 1108–1114, DOI: [10.1002/sia.6284](https://doi.org/10.1002/sia.6284) (2017).
34. Rondin, L. *et al.* Nanoscale magnetic field mapping with a single spin scanning probe magnetometer. *Appl. Phys. Lett.* **100**, 153118, DOI: [10.1063/1.3703128](https://doi.org/10.1063/1.3703128) (2012).
35. Wu, Y. *et al.* Molecular rheometry: direct determination of viscosity in Lo and Ld lipid phases via fluorescence lifetime imaging. *Phys. Chem. Chem. Phys.* **15**, 14986, DOI: [10.1039/c3cp51953h](https://doi.org/10.1039/c3cp51953h) (2013).
36. Drexhage, K. Influence of a dielectric interface on fluorescence decay time. *J. Lumin.* **1-2**, 693–701, DOI: [10.1016/0022-2313\(70\)90082-7](https://doi.org/10.1016/0022-2313(70)90082-7) (1970).
37. Costa, L. *et al.* Spectroscopic Investigation of Local Mechanical Impedance of Living Cells. *PLoS ONE* **9**, e101687, DOI: [10.1371/journal.pone.0101687](https://doi.org/10.1371/journal.pone.0101687) (2014).
38. Rico, F., Gonzalez, L., Casuso, I., Puig-Vidal, M. & Scheuring, S. High-Speed Force Spectroscopy Unfolds Titin at the Velocity of Molecular Dynamics Simulations. *Science* **342**, 741–743, DOI: [10.1126/science.1239764](https://doi.org/10.1126/science.1239764) (2013).
39. Sader, J. E. *et al.* A virtual instrument to standardise the calibration of atomic force microscope cantilevers. *Rev. Sci. Instruments* **87**, 093711, DOI: [10.1063/1.4962866](https://doi.org/10.1063/1.4962866) (2016).
40. Proksch, R., Schäffer, T. E., Cleveland, J. P., Callahan, R. C. & Viani, M. B. Finite optical spot size and position corrections in thermal spring constant calibration. *Nanotechnology* **15**, 1344–1350, DOI: [10.1088/0957-4484/15/9/039](https://doi.org/10.1088/0957-4484/15/9/039) (2004).
41. Milhiet, P. E., Vié, V., Giocondi, M.-C. & Le Grimallec, C. AFM Characterization of Model Rafts in Supported Bilayers. *Single Mol.* **2**, 109–112, DOI: [10.1002/1438-5171\(200107\)2:2<109::AID-SIMO109>3.0.CO;2-L](https://doi.org/10.1002/1438-5171(200107)2:2<109::AID-SIMO109>3.0.CO;2-L) (2001).
42. Redondo-Morata, L., Giannotti, M. I. & Sanz, F. Influence of Cholesterol on the Phase Transition of Lipid Bilayers: A Temperature-Controlled Force Spectroscopy Study. *Langmuir* **28**, 12851–12860, DOI: [10.1021/la302620t](https://doi.org/10.1021/la302620t) (2012).
43. Gumí-Audenis, B. *et al.* In-plane molecular organization of hydrated single lipid bilayers: Dppc: cholesterol. *Nanoscale* **10**, 87–92, DOI: [10.1039/C7NR07510C](https://doi.org/10.1039/C7NR07510C) (2018).