

A study of the indirect geometry cold neutron spectrometer, BIFROST, from neutron source position to detector position

M. Klausz^{a,b,c,*}, K. Kanaki^b, T. Kittelmann^b, R. Toft-Petersen^{b,d}, J.O. Birk^e,
M.A. Olsen^e, P. Zagvyai^{a,c}, R.J. Hall-Wilton^{b,f}

^a*Hungarian Academy of Sciences, Centre for Energy Research, 1525 Budapest 114., P.O. Box 49., Hungary*

^b*European Spallation Source ESS ERIC, P.O Box 176, SE-221 00 Lund, Sweden*

^c*Budapest University of Technology and Economics, Institute of Nuclear Techniques, 1111 Budapest, Műegyetem rakpart 9., Hungary*

^d*Technical University of Denmark, Department of Physics, DK-2800 Kongens Lyngby, Denmark*

^e*Nanoscience Center, Niels Bohr Institute, University of Copenhagen, DK-2100 Copenhagen Ø, Denmark*

^f*Università degli Studi di Milano-Bicocca, Piazza della Scienza 3, 20126 Milano, Italy*

Abstract

The European Spallation Source (ESS) is going to be the most powerful spallation neutron source in the world and the flagship of neutron science in the upcoming decades. The exceptionally high neutron flux will provide unique opportunities for scientific experiments, but also set high requirements for the detectors. One of the most challenging aspects is the rate capability and in particular the peak instantaneous rate capability, i.e. the number of neutrons hitting the detector per channel or cm² at the peak of the neutron pulse.

The primary purpose of this paper is to estimate the incident rates that are anticipated for the BIFROST instrument planned for ESS, and also to demonstrate the use of powerful simulation tools for the correct interpretation of neutron transport in crystalline materials.

A full simulation model of the instrument from source to detector position, implemented with the use of multiple simulation software packages is presented. This paper includes the first application of a new pyrolytic graphite model

*Corresponding author

Email address: milan.klausz@energia.mta.hu (M. Klausz)

and a comparison of different simulation tools to highlight their strengths and weaknesses.

Keywords: Geant4, McStas, neutron detector, neutron scattering, neutron spectroscopy, crystallography

1. Introduction

The European Spallation Source (ESS) ERIC [1, 2, 3] is designed to operate using the most powerful spallation neutron source in the world, and to provide unprecedentedly high neutron fluxes for instruments of various neutron techniques. One of these instruments is BIFROST [4, 5], a high flux, indirect geometry, cold spectrometer, optimised for small samples and extreme environments. BIFROST provides exceptional flux as it can operate in a white beam mode. This flux allows for entirely new options for detailed investigations of complex multimode dynamics, hybrid modes, electro-magnons, spin wave continua and gap studies, under extreme conditions with controlled temperature, pressure, and magnetic fields.

Harnessing the full ESS pulse by employing a polychromatic beam carries a high potential risk for detector rates that can saturate the detectors and therefore degrade the performance of the instrument. Position sensitive ^3He tubes are the “gold standard” for neutron detection [6]. They are however quite rate limited. Non-position sensitive tubes saturate at 100 kHz; however for position sensitive ^3He tubes, operation at instantaneous rate above 30 kHz can be problematic. The exact rate capability of a detector is dependent on the details of readout electronics. It is therefore essential to evaluate the rates anticipated for high-flux instruments [7, 8, 9], in order to extract the respective detector requirements.

Monte Carlo simulation plays a key role in the development and characterisation of instruments as a reliable, cheap and versatile tool [10]. Feedback from simulations taken into account in the development of the instrument design can reduce the number of physical prototypes needed, and also enables the quan-

tification of otherwise unmeasurable properties. This is particularly the case for complicated instruments, such as BIFROST. Development of complete and detailed instrument simulation models enables simulations from source to detectors, offering the opportunity to discover and decouple otherwise undetectable cumulative effects. These models can provide valuable input for developing calibration and correction routines for data reduction and analysis, and could later be used for experiment planning by users, to predict experimental conditions from specific proposed samples, i.e. sample size and composition.

To make such full instrument simulations possible [7, 11], advanced simulation tools like NCrystal [10, 12] have been developed, which enables Monte Carlo simulations of thermal neutrons in crystals, and Monte Carlo Particle List (MCPL) [13, 14], which enables communication between different software packages. These tools can greatly enhance the capabilities of the existing and widely used simulation software like McStas [15, 16] and Geant4 [17, 18, 19].

In this study multiple Monte Carlo simulation tools are used together to implement a full simulation model of the BIFROST instrument from the neutron source to the detector position. This full model is then used to estimate the incident detector rates that are anticipated in the case of the highest possible incident neutron intensity – that will be consequently mentioned as ‘worst-case’ – and normal-use scenarios. Based on these studies the detector requirements for rate capability can be extracted.

In the following sections, the instrument and simulation models geometries and tools are introduced first, followed by the presentation of incident rates for elastic peaks for various instrument parameters, sample types, sizes and mosaicities, along with a demonstration of the differences between McStas and Geant4 simulation results. The study concludes with the demonstration of the elastic signal of a standard calibration sample.

2. BIFROST instrument and simulation model

2.1. The BIFROST instrument

BIFROST is a 162 m long cold neutron spectrometer intended to be built as a first tranche instrument for ESS. It combines an indirect geometry time-of-flight (ToF) front end, and an angular and energy multiplexed crystal analyser-based back end. A back end similar in principle to that installed recently at the CAMEA spectrometer at the PSI [20, 21]. BIFROST is designed [22] to maximise the use of the ESS long pulse to enable measurements on small samples and study dynamic properties, under extreme conditions with controlled temperature, pressure, and magnetic fields. The envisioned application fields include materials science, magnetism, life sciences and planetary sciences [5].

The instrument consists of three main technical subsystems: the beam transport and conditioning system, the sample exposure system and the scattering characterisation system. The schematic model of the instrument is depicted in Fig. 1.

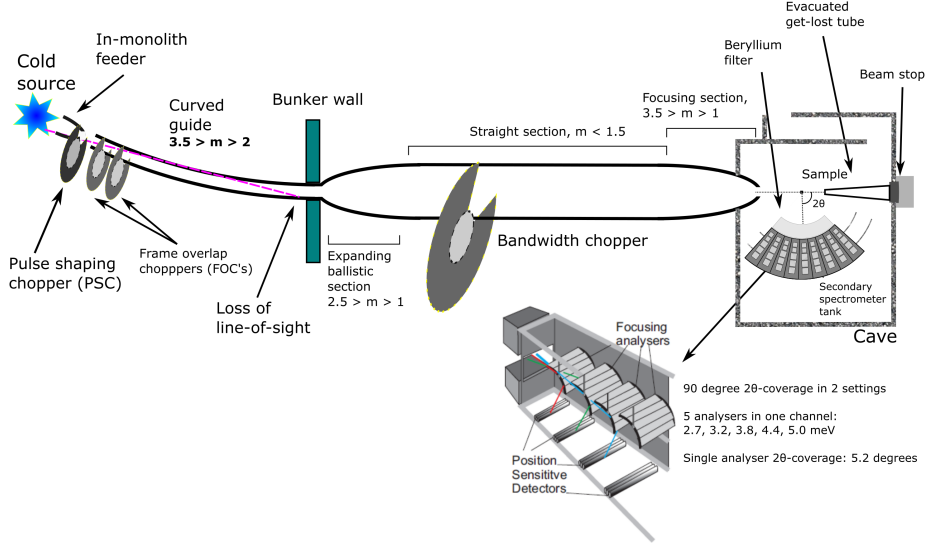


Figure 1: Schematic model of BIFROST from source to detection position. From [23].

The beam transport and conditioning system is relatively simple. It has a

curved guide section inside the bunker to lose line-of-sight, and four choppers as the only moving parts. Three of these choppers are placed inside the bunker and the fourth one is placed in the middle of the instrument. The first one is a pulse-shaping chopper, that is the only chopper determining the energy resolution. The other three – two Frame Overlap and one Bandwidth Chopper – serve to sort out unwanted frames from the fast pulse-shaping chopper and to avoid pulse overlap at the sample position respectively. The pulse-shaping chopper can reduce the ESS pulse width by a factor of up to 30 to match the best analyser resolution, or allow the full pulse to reach the sample that will result in a relaxed resolution but an order of 10^{10} n/s/cm² flux on sample. It is this mode which poses the greatest rate challenge for the detectors.

The sample exposure system allows measurements with strong magnetic field, high pressure and cryogenic temperatures. One of the main limitations today in single crystal neutron spectroscopy is that measurements are only possible with large samples, which are not available for many sample types, but BIFROST will enable the study of sub-cubic millimetre samples thanks to its exceptional flux on sample and the efficient scattering characterisation system.

The scattering characterisation system in Fig. 2 consists of the filtering system and the secondary spectrometer tank, covering a 90° scattering angle in the horizontal plane in two tank settings. The filtering system, that is essential for background reduction on BIFROST, includes a cooled beryllium filter with roughly 90% transmission of neutrons with an energy below 5 meV (4.05 Å) but very low transmission of neutrons with energies above, and coated lamellas as a radial collimator [24].

The secondary spectrometer tank houses multiple sets of analysers and detectors for different neutron energies, arranged in nine so called Q-channels. Depending on the scattering angle of a neutron on the sample, it enters one of the Q-channels, which are separated by cross-talk shielding between them. In each Q-channel, several crystal analyser arcs are placed one behind the other to select different final energies by scattering neutrons vertically (down) towards the corresponding set of position sensitive detectors, employing Rowland focus-

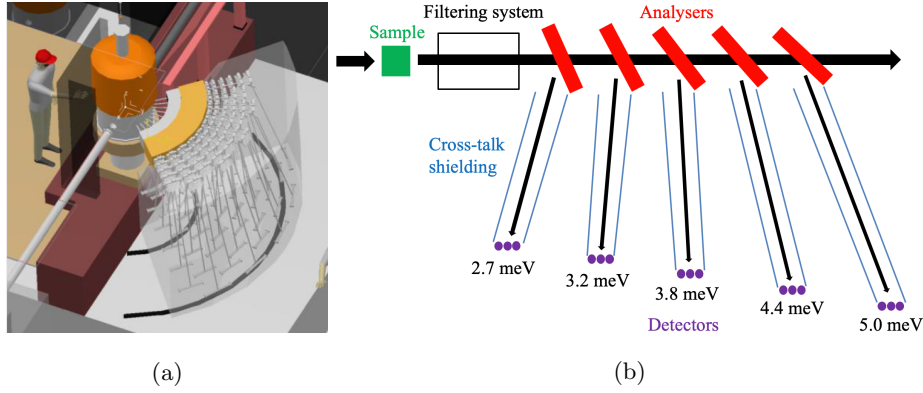


Figure 2: The BIFROST scattering characterisation system. 3D model of the secondary spectrometer tank from [4] (a) and side view sketch of the sample and a single Q-channel with five energy channels in it (b).

ing (see Fig. 2b). Further cross-talk shielding is applied to ensure that neutrons can reach the detectors only by scattering from the corresponding analyser arcs.

With this arrangement, BIFROST utilises a variant of a novel analyser setup called CAMEA [21, 20], an acronym for continuous angle multiple energy analysis. Enabling multi-energy analysis in a single Q-channel by placing the analysers for higher neutron energies behind the ones for lower energies is possible due to the high transparency of the 1 mm thin highly-oriented pyrolytic graphite blades [25]. The blades to be used have high mosaicity (60 arcmin) to apply the prismatic analyser concept [26], using ^3He detector triplets for all five neutron energies chosen for BIFROST (2.7 meV, 3.2 meV, 3.8 meV, 4.4 meV, 5.0 meV) in each Q-channel. According to the prismatic analyser concept, each of the three detectors of a triplet records a slightly different region of energy due to the high analyser mosaicity, as neutrons with different energies are scattered in slightly different directions.

In order to provide enough space for the detector tubes, the analysers and corresponding detectors in adjacent Q-channels are slightly shifted radially. The sample-analyser distances are increased or decreased by 4.6–7.5% in two out of three Q-channels, however the analyser-detector distances are kept unchanged to

keep the detectors of same energy on the same vertical planes and by that keeping the spectrometer tank geometry simple. As a result, the sample-analyser distance is shorter or longer than the analyser-detector distance in two out of three Q-channels, showing a slight asymmetry to the Rowland-geometry. The three different type of Q-channels are repeated three times, giving the nine Q-channels a “triple stagger” geometry.

The simulation of the BIFROST instrument is divided into two parts: the simulation of the 162 m long beam transport and conditioning system until the end of the last guide section, and the simulation of the sample and scattering characterisation system together (see Fig. 3). The first part is done using McStas only; the second part is implemented and simulated in both McStas and Geant4, in order to compare the results of these simulation tools and to demonstrate why it is advantageous to use Geant4 for the back-end of the instrument. The transition between the two parts is facilitated by the aforementioned MCPL tool, that is a binary file format dedicated for storage and interchange of particles between various Monte Carlo simulation applications, like McStas, Geant4, McXtrace [27] and MCNP [28]. For the simulation of neutron transport in crystalline materials, NCrystal is used in both McStas and Geant4. In the next subsections all simulation tools and models are introduced.

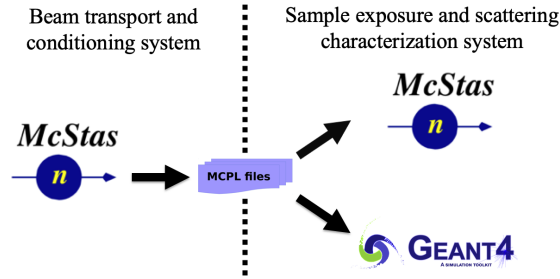


Figure 3: Outline of simulation scheme.

2.2. McStas model

McStas [15, 16] is a Monte Carlo simulation tool dedicated for simulation of neutron scattering instruments and experiments. It is user-friendly, cross-

platform, open source and uses a ray-tracing algorithm which enables fast neutron transport simulations over long distances and through many components, that is necessary for long instruments like BIFROST.

For the simulation of neutrons from the source to the end of the beam transport and conditioning system, a previously developed available McStas model of the instrument is used [29]. This model, depicted in Fig. 4 contains the butterfly moderator source (“ESS_butterfly”) [30, 31], the four choppers, all guide sections, and several McStas monitor components to characterise the beam at multiple locations along the guide. The source is used with the highest intensity, so the deduced rate numbers in this paper correspond to the maximum accelerator power of 5 MW. Expected rates will scale linearly with source power for constant proton energy. At the end of the last guide section all neutron data are saved in an MCPL file using the “MCPL_output” McStas component. This file serves as input for both the McStas and the Geant4 simulation model of the second part of the instrument.

The McStas model of the sample and scattering characterisation system, depicted in Fig. 4, contains a crystalline sample, one Q-channel including all five analyser arcs, and McStas monitor components at several places probing ToF, energy and position distribution of neutrons, in order to examine the change of the neutron beam. The analyser arcs consist of 7–9 blades using the “NCrystal_sample” component with pyrolytic graphite material, described in more detail in section 2.4. NCrystal is also used for all crystalline samples throughout this study. The simulation model does not contain the sample environment, eight out of the nine Q-channels, the filtering system, cross-talk shielding and the detectors.

Using a reduced geometry and excluding any model of the sample environment is intentional, aiming to get a conservative estimate in terms of the highest detector rates, but implementing only one Q-channel is the result of a limitation coming from the linearity of the McStas simulation process. In a McStas instrument definition file, the geometrical components like the source, guide sections, choppers, slits and the sample are placed one after the other. McStas by default

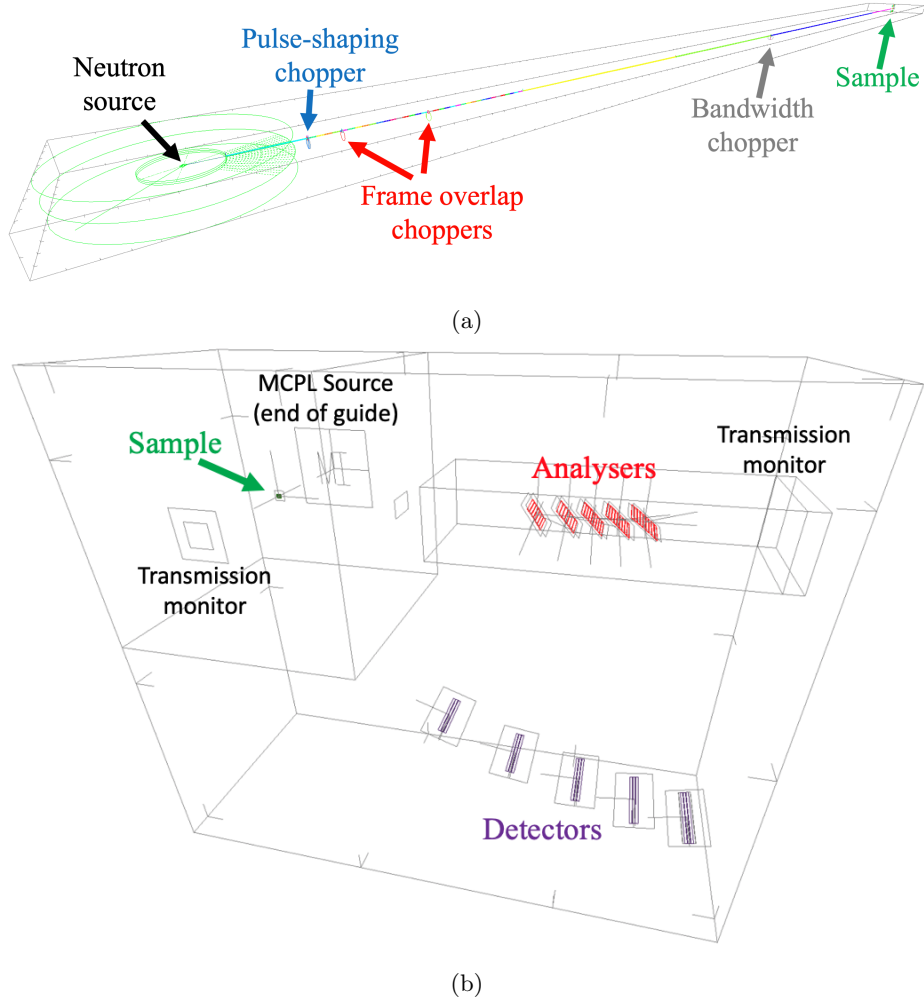


Figure 4: McStas model of the beam transport and conditioning system (a) and the sample and scattering characterisation system (b). The figures are at different scale.

propagates neutrons from component to component in the exact order as they appear in this file. All neutrons that miss or do not interact with the component downstream are removed from further simulation. This process makes the simulation of long instruments really fast, but on the other hand restricts the neutrons to follow one exact path, that does not allow the simulation of multiple Q-channels simultaneously. It is possible to change this behaviour by grouping

components together, as described in the users and programmers guide [32]. This way it is possible to some extent to split the beam, by having a group of components as the potential next target of neutrons, but after the interaction with one of the group members, McStas tries to propagate all neutrons to the same component that appears next in the instrument file. This means that if a user wants to split the beam and allow propagation in multiple directions through different components, consecutive groups must be implemented, all of which include the subsequent component in each direction. Extensive use of such grouping makes the instrument file immensely complex and still prohibits multiple interactions within one group, or back and forth propagation between groups of components.

For these reasons, only one Q-channel is implemented within which this technique is used to handle the five sets (arcs) of analyser blades, all of which divide the beam into the partition that is scattered towards the detectors and the partition that is propagated towards the next set of analysers (or the beam stop behind the last set). In order to allow neutrons to proceed without interaction with a set of analyser blades – due to transmission or simply avoiding them all – an extra virtual component is added to each group, that mimics an interaction without changing the neutron state, and thereby prevents neutrons to be removed from the simulation. As mentioned, neutrons still have to follow the order of the groups, so back-scattering or multiple scattering among blades of the same arc is still not possible in the simulations.

Although the cross-talk shielding between energy- and Q-channels is not explicitly included in the model, as a consequence of the above described process, a neutron can reach a particular detector tube only by scattering in one of the corresponding analyser blades. This is practically equivalent to an ideal cross-talk shielding absorbing all stray neutrons. The case is similar for the filtering system, that is replaced by a monitor component, that transmits all neutrons below 7 meV energy and none above. As mentioned earlier, the transmission of beryllium drops sharply around 5 meV, that is in fact the highest of the five final energies selected by the analysers. Simulation of effects of this transition in the

transmission is out of the scope, and using ideal transmission in the 0–7 meV energy region keeps the rate estimates conservative.

The intent is to determine the incident detector rates, therefore the simulation of the detection process is out of scope. Detectors are modelled with McStas monitor components, and a neutron is counted as incident for a detector tube if it crosses the plane at the centre of the detectors within the outline of that particular tube. The sample–analyser distance is equal to the analyser–detector distance, meaning that symmetrical Q-channels are modelled.

As this subsection demonstrates, using McStas to model such complex system as the analyser-detector system of BIFROST in detail is not straightforward and is subject to certain limitations.

McStas version 2.5 is used for the simulations.

2.3. Geant4 model

Geant4 [17, 18, 19] is a general purpose Monte Carlo particle transport toolkit developed at CERN with applications in many fields, e.g. high energy physics, nuclear physics, accelerators and medical physics. Its usability for simulation of neutron detectors has been greatly improved by the ESS Detector Group by building a framework [10, 33] around it, which adds several functionalities and integrates NCrystal and MCPL.

The Geant4 simulation model of the sample and scattering characterisation system, depicted in Fig. 5, contains the same parts as the McStas model (crystalline sample, all five analyser arcs in a Q-channel) but with the option to simulate with all Q-channels included. To make the results comparable with the McStas model, the McStas monitor components are mimicked with empty volumes with exactly the same location and surface, in order to create histograms with the same predefined spatial, energy and time-of-flight (ToF) resolution.

Although the cross-talk shielding could be easily implemented in the model, for the same comparison purposes it is replaced by certain conditions applied at data in the analysis level. This means that neutrons cannot skip parts of the model. They can, however, scatter back and forth between the geometrical

components many times, since there is no such limitation on order of length for the path of the neutrons in Geant4, as there is in McStas. This gives the possibility to analyse effects of cross-talk on signals and to evaluate shielding strategies.

The sample–analyser distances in the Q-channels of the simulation model correspond to the real geometrical parameters. The analyser–detector distances, however, are the same as the sample–analyser distances in all Q-channels, meaning that the slight asymmetry of the non-symmetrical Q-channels described in section 2.1 is neglected. This only affects section 5 where all nine Q-channels are used, in all other simulations a single symmetrical Q-channel is modelled.

The Geant4 physics list used is QGSP_BIC_HP_EMZ that uses high precision models and cross sections for neutron energies lower than 20 MeV, and allows the correct treatment of thermal and cold neutrons when combined with NCrystal.

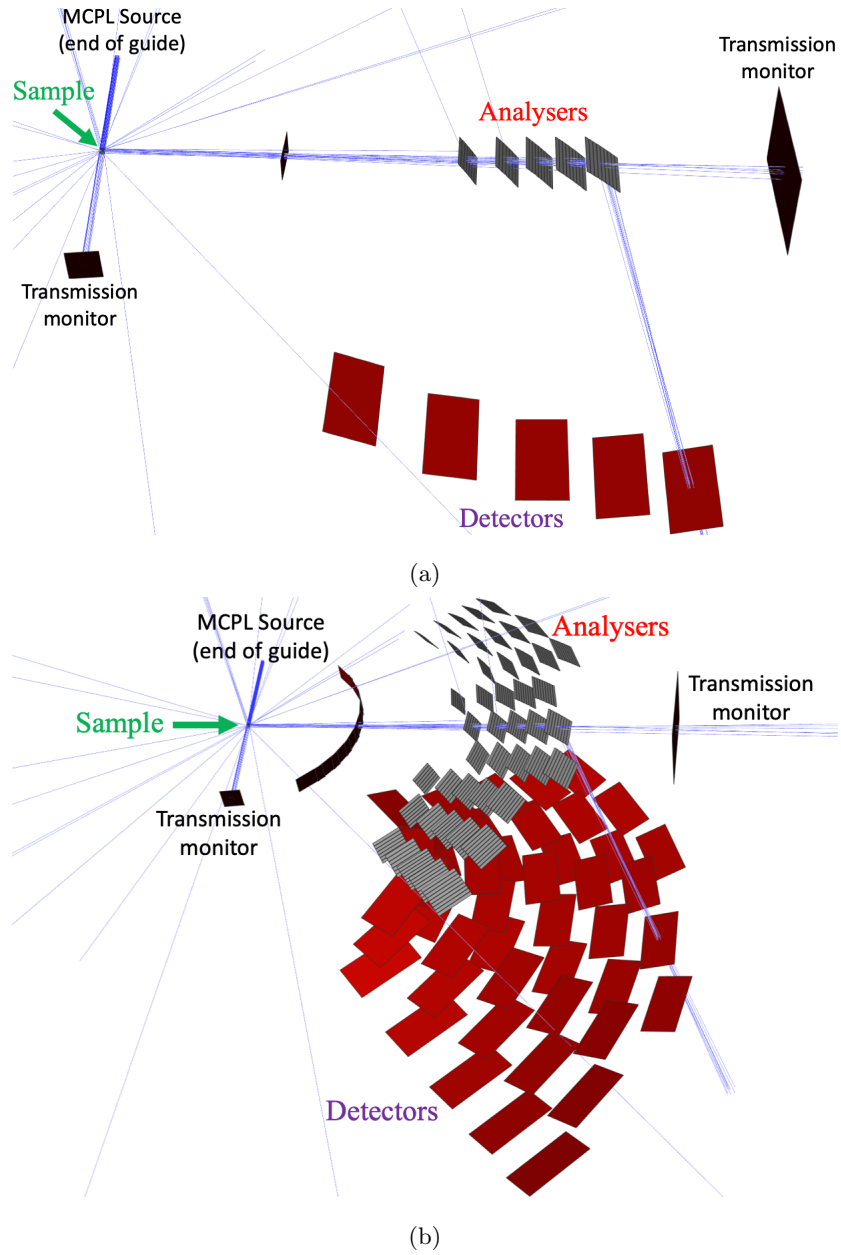


Figure 5: Geant4 model of the sample and scattering characterisation system including one Q-channel (a) or all Q-channels (b). The blue lines indicate simulated neutron paths.

2.4. *NCrystal*

NCrystal [10, 12] is a novel open source software package for modelling thermal neutron transport in crystalline materials. It consists of a data library and associated tools which enables calculations for Monte Carlo simulations. It can be used together with McStas and Geant4, to enhance their capabilities of the correct treatment of neutron transport in typical components of neutron instruments, including beam filters, monochromators, analysers, samples and detectors. Physics modelled by NCrystal includes both coherent elastic (Bragg diffraction) and incoherent or inelastic (phonon) scattering. It treats all valid Bragg diffractions on each reflection plane explicitly and is able to use various models for inelastic scattering on phonons. Its data library already contains the most popular crystals and the results are validated against the EXFOR database and existing crystallographic software.

NCrystal focuses initially on scattering in single-crystals or polycrystalline materials and powders. Most single crystalline materials are appropriately modelled with crystallites orientated around some reference orientation with a Gaussian distribution that has a standard deviation of the mosaicity of the crystal. There are, however, single crystalline materials with crystallite distributions so different from Gaussian that this approximation does not hold. One of these materials is pyrolytic graphite, that is widely used as a monochromator and analyser in neutron instruments. This is precisely the case for BIFROST, where 369 highly oriented pyrolytic graphite analyser blades are used altogether in the nine Q-channels. Graphite has a layered structure, made up of graphene sheets in which carbon atoms are arranged in a hexagonal lattice. In highly oriented pyrolytic graphite the crystallite axes orthogonal to the graphene sheets are distributed along a preferred direction, suitable for description with a Gaussian mosaicity distribution, but the orientation around this axis is completely random, resulting in powder-like features in neutron scattering. Recent developments enabled NCrystal to handle materials with this kind of structure by using a specialised model for layered crystals.

As mentioned earlier, NCrystal is used for the crystalline materials both in

the McStas and Geant4 model of BIFROST. The samples are modelled as single-crystals, but for the pyrolytic graphite analysers the new layered crystal model is used. Demonstrating some of the NCrystal pyrolytic graphite properties, the components of total cross-sections and a distribution of randomly sampled scatter angles are depicted in Fig. 6 for powder sample, not the layered crystal distribution.

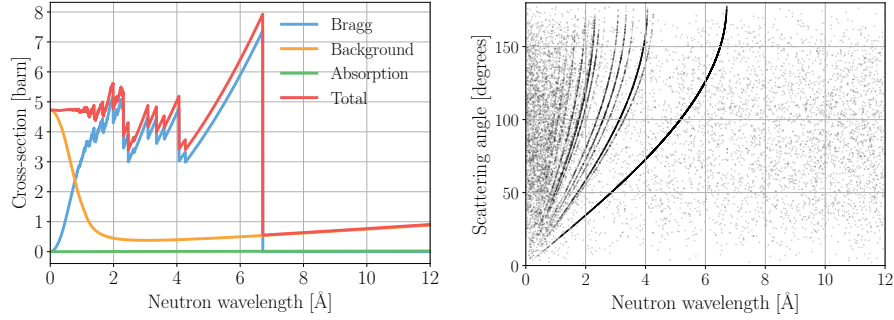


Figure 6: Components of the total interaction cross-section (left) and two-dimensional scatter-plot (right) for NCrystal pyrolytic graphite powder.

Up until the latest version (v1.0.0), NCrystal treats absorption with the simple model of absorption cross-sections being inversely proportional to the neutron velocity. The absorption cross-section for a particular neutron velocity is calculated by scaling the value given at the reference velocity of 2200 m/s. This applies for NCrystal materials used in McStas simulations, but not for those in Geant4 simulations. Geant4 models secondary particles produced in absorption, therefore the NCrystal plugin does not interfere with the Geant4 absorption physics at all. As an example, the minor differences in absorption cross-section for pyrolytic graphite is shown in Fig. 7.

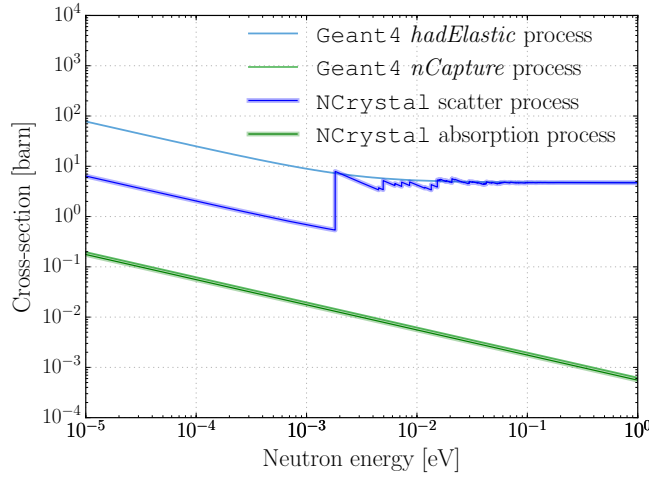


Figure 7: Comparison of neutron interaction cross-sections of pyrolytic graphite powder in Geant4 and NCrystal. Lines are barely distinguishable for Geant4 neutron capture and NCrystal absorption cross-sections, the former being higher in the whole energy range depicted.

All results of simulations have systematic and statistical uncertainties. Unlike systematic uncertainties which are present mainly due to imperfect modeling of the system, statistical uncertainties can be reduced by increasing the number of simulated particles. Throughout this work where the uncertainty is not significant, it is not indicated. Sources of systematic uncertainties are considered and not expected to change the conclusions.

With the simulation models and tools described in this section, a detailed analysis is carried out, in order to give an estimate of the incident rates that are anticipated for detectors at BIFROST. The incident detector rates for elastic peaks using various samples and instrument parameters are presented in the next section, as well as a comparison of results with McStas and Geant4 simulations.

3. Incident rates for coherent elastic peaks

Determination of anticipated detector rates for an instrument is a key part of defining requirements for the detectors to be used. It can prevent the detector rate capability from becoming the bottleneck of experiments or a source of performance degradation. For this reason, the simulation tools and models described in the previous section are used to determine the highest time-averaged and the highest instantaneous (peak) incident rates for the detector tubes. These rates are determined for the worst-case and for more realistic operational conditions.

The worst-case incident rate for a single detector tube occurs when a strong Bragg peak from a single crystal sample gets reflected to it. To get such elastic scattering on the sample, Bragg's law must apply:

$$2 \cdot d_{hkl} \sin(\theta) = n \cdot \lambda, \quad (1)$$

where d_{hkl} is the interplanar distance of the crystal, θ is the scattering angle, n is a positive integer and λ is the wavelength of the incident neutrons. To get neutrons scattered on the sample onto the detectors, their energy has to match one of the energies selected by the analyser arcs. Based on the ESS source spectrum, the guide transmission and the energy resolution of the analysers, the highest incident rates are expected for the 5.0 meV (4.045 Å) neutrons. This sets the requirement of the lattice parameter being larger than 4.045/2 Å, according to Eq. 2, derived from Eq. 1.

$$d_{hkl} \geq \frac{\lambda}{2}. \quad (2)$$

Regarding highly reflective materials that would result in the highest detector rates, the truly worst-case sample would be a pyrolytic graphite single crystal ($d_{002}=3.3555$ Å), but to get results from a less unrealistic sample with large enough lattice parameter and strong Bragg peak, simulations are also done using an yttrium oxide (Y_2O_3) single crystal ($d_{2-2-2}=3.0724$ Å).

Further parameters that influence the rates on detectors are the pulse-shaping chopper opening time, sample size, and sample mosaicity. To get the highest possible rates, maximum flux mode is used, when the pulse-shaping chopper is fully open, resulting to a 10^{10} n/s/cm² flux on sample. The instrument is designed to facilitate measurements on small samples, but dimensions up to 1.5 cm are possible, therefore cylindrical samples with the diameter and height of 1.5 cm are used in the simulations. The highest rates are expected when the mosaicity of the sample matches that of the analysers, so the sample mosaicity is set to 60 arcmin.

In order to realise the simulations, the samples are oriented to fulfil the Bragg law criteria for the incoming 5 meV neutrons on the selected scattering planes, and the single Q-channel modelled is rotated according to the resulting scattering angle. In BIFROST the whole scattering characterisation system can be rotated around the sample, so having a Bragg peak from any sample in the exact direction of a single Q-channel is perfectly realistic.

The simulation with the parameters described above is done in two steps. First, the neutron transport from the source to the end of the beam transport and conditioning system is done with the McStas model, saving neutron data at the end in an MCPL file. This file is then used as the source term for the simulation of the sample and scattering characterisation system, that is done with both McStas and Geant4. The results are presented in the following subsection (3.1). In the subsequent subsection the impact of different parameters like sample and analyser mosaicity, sample size, pulse-shaping chopper opening time is scanned, but due to the good agreement of McStas and Geant4 simulations results (demonstrated later in subsection 3.1), this is only done using the Geant4 model.

3.1. McStas – Geant4 comparison

3.1.1. Pyrolytic graphite sample

Fig 8 depicts the time-averaged energy spectra of neutrons at the sample and different parts of the scattering characterisation system from the McStas and

Geant4 simulations using a pyrolytic graphite single crystal sample. Due to the large mosaicity of the sample, neutrons of a broad energy range – centred around 5.0 meV – are scattered toward the single Q-channel, reaching the 2.7 meV analysers and therefore missing from the sample transmission spectrum. The change of the spectrum between the 2.7 meV and 5.0 meV analysers is caused by the spread of neutrons, and absorption in the analysers. The neutrons selected by the 5.0 meV analysers are scattered towards the corresponding detector triplet, and are therefore absent from the analyser transmission spectrum.

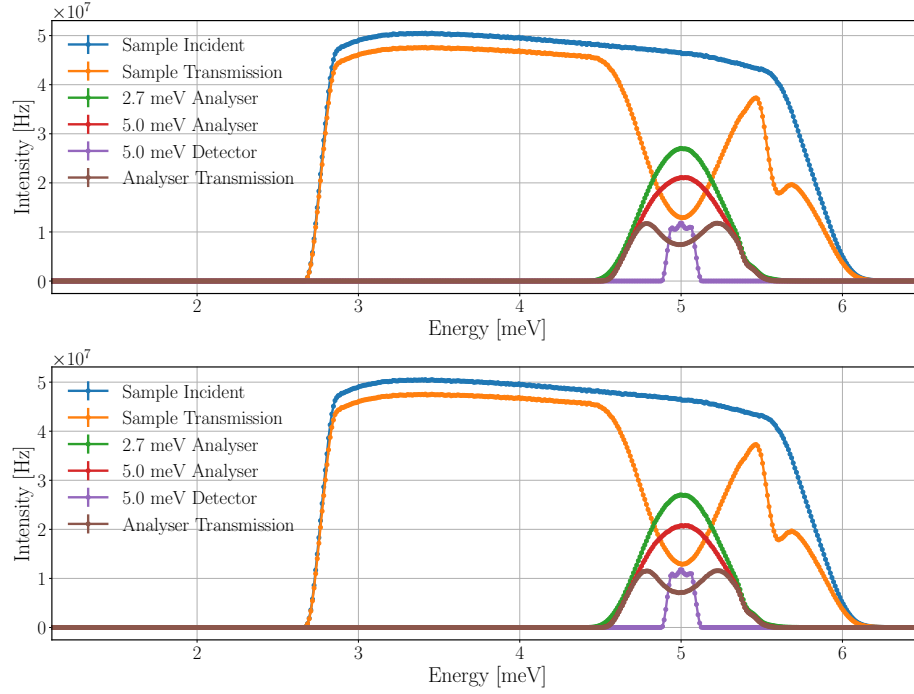


Figure 8: Time-averaged neutron energy spectra at the sample and the scattering characterisation system with pyrolytic graphite single crystal sample in McStas (top) and Geant4 (bottom) simulation. Incident beam on sample (in blue), beam transmitted through the sample (in orange), beam on the set of analysers for 2.7 meV neutrons (in green), beam on the set of analysers for 5.0 meV neutrons (in red), neutrons hitting the detector triplet for 5.0 meV (in purple), beam transmitted through all sets of analysers (in brown). The lines are only joining the points.

The time-averaged neutron intensities acquired by the integration of the

energy spectra are shown in Tab. 1. The results of the McStas and Geant4 simulations agree with only minor differences.

| Position | McStas [Hz] | Geant4 [Hz] |
|-----------------------|----------------------|----------------------|
| Sample Incident | $1.45 \cdot 10^{10}$ | $1.45 \cdot 10^{10}$ |
| Sample Transmission | $1.16 \cdot 10^{10}$ | $1.16 \cdot 10^{10}$ |
| 2.7 meV Analyser | $1.50 \cdot 10^9$ | $1.50 \cdot 10^9$ |
| 5.0 meV Analyser | $1.20 \cdot 10^9$ | $1.19 \cdot 10^9$ |
| 5.0 meV Detector | $1.98 \cdot 10^8$ | $1.98 \cdot 10^8$ |
| Analyser Transmission | $7.47 \cdot 10^8$ | $7.31 \cdot 10^8$ |

Table 1: Time-averaged neutron intensities at the sample and different parts of the scattering characterisation system with pyrolytic graphite single crystal sample with McStas and Geant4.

The structure of the energy spectrum of the detector triplet is the result of summing the spectra of all three detectors, as depicted in Fig 9. As mentioned earlier, the analysers scatter neutrons with slightly different energies in different directions due to their high mosaicity, according to the prismatic analyser concept. As a result of this vertical spread, the three detectors of the triplet record slightly different regions of energy. It is not only the analysers, that spread the beam due to their high mosaicity, but the sample too. The sample spreads the neutron beam horizontally, by scattering neutrons with higher energies in larger, and neutrons with lower energies in smaller angle than the scattering angle given by the Bragg formula (Eq. 1). The combined effect of these processes is visible in Fig 10, showing a diagonal shape in the time-averaged neutron intensities in the plane over the detector tubes.

Integrating the incident neutron intensities in Fig. 10 over the areas of the tubes gives the time-averaged incident rates for the tubes, presented in Tab. 2. The results show that the time-averaged incident rate for a single detector tube can be almost as high as 70 MHz.

Given that BIFROST is a ToF instrument at a spallation neutron source, the incident detector rate naturally has a pulsed time structure. The ToF dis-

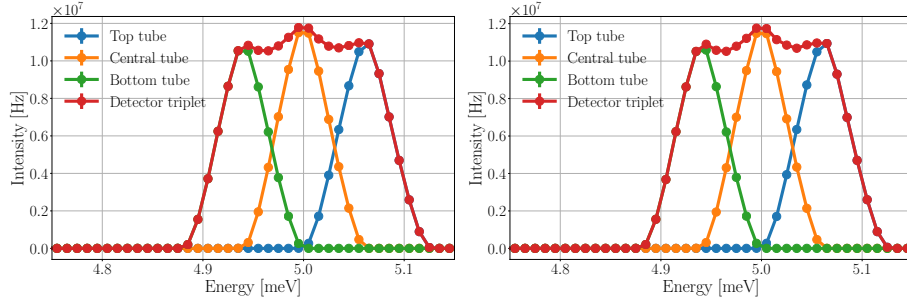


Figure 9: Time-averaged incident neutron energy spectra of the 5 meV detector tubes in McStas (left) and Geant4 (right) simulations, with pyrolytic graphite sample. The lines are only joining the points.

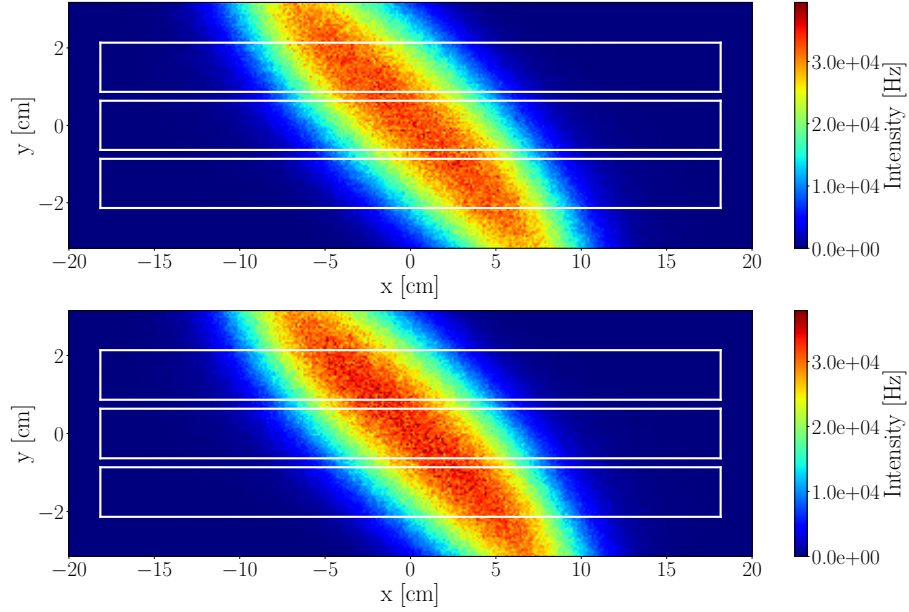


Figure 10: Time-averaged neutron intensities in the plane of the 5 meV detector tubes in McStas (at top) and Geant4 (at bottom) simulations, with pyrolytic graphite sample. The white lines indicate the outlines of the detector tubes.

tribution of a single pulse on the 5 meV detectors is depicted in Fig. 11. By taking into account only those neutrons which arrive to the detectors at the peak of their ToF distribution in a short time range of 0.1 ms (see Fig. 12), the integration over the areas of the detector tubes gives the peak incident rates,

| Detector tube | McStas [Hz] | Geant4 [Hz] |
|---------------|-------------------|-------------------|
| Top | $6.69 \cdot 10^7$ | $6.70 \cdot 10^7$ |
| Central | $6.95 \cdot 10^7$ | $6.93 \cdot 10^7$ |
| Bottom | $6.20 \cdot 10^7$ | $6.20 \cdot 10^7$ |

Table 2: Time-averaged incident neutron rates of the 5 meV detector tubes in McStas and Geant4 simulations, with pyrolytic graphite sample.

as presented in Tab. 3. Due to the distinct energy range and therefore different ToF spectra of the tubes, the highest peak incident rate occurs at different times for each tube of a triplet. The results demonstrate that the peak incident rate on a single detector tube can be as high as 1.7 GHz.

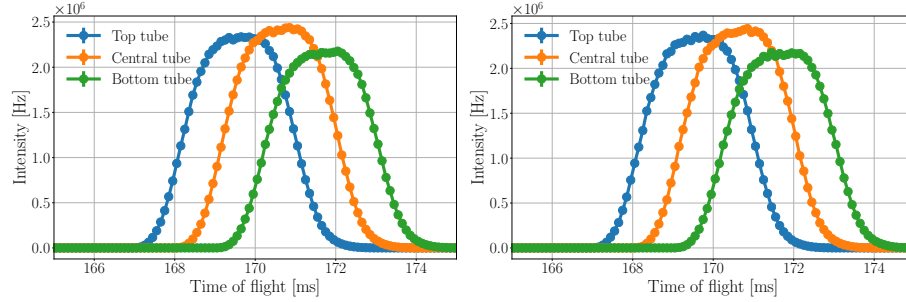


Figure 11: ToF spectrum of neutrons at the 5 meV detector triplet in McStas (left) and Geant4 (right) simulations, with pyrolytic graphite sample. The lines are only joining the points.

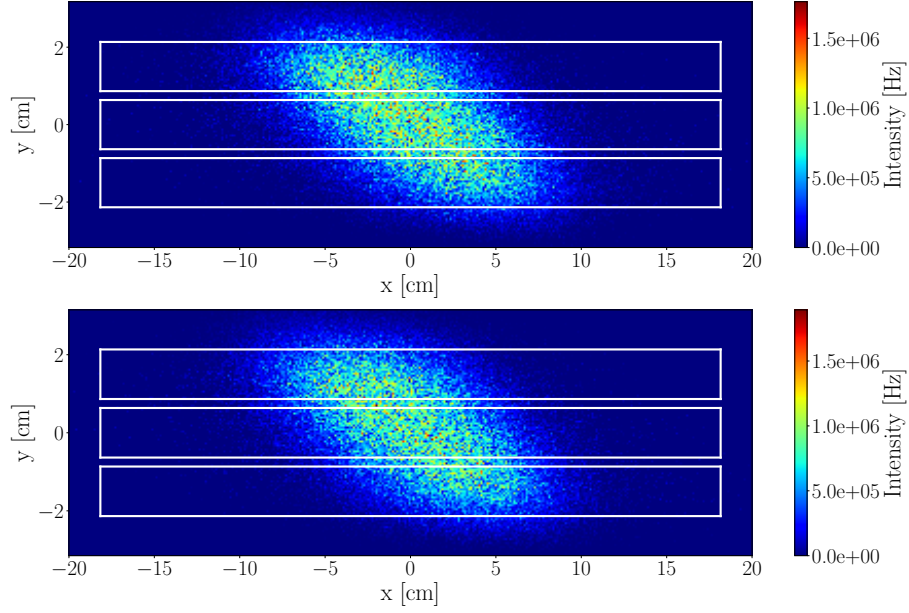


Figure 12: Neutron intensities in the plane of the 5 meV detector tubes at the peak of the ToF spectrum of the central tube in McStas (top) and Geant4 (bottom) simulations, with pyrolytic graphite sample. The white lines indicate the outlines of the detector tubes.

| Detector tube | McStas [Hz] | Geant4 [Hz] |
|---------------|-------------------|-------------------|
| Top | $1.67 \cdot 10^9$ | $1.69 \cdot 10^9$ |
| Central | $1.74 \cdot 10^9$ | $1.75 \cdot 10^9$ |
| Bottom | $1.56 \cdot 10^9$ | $1.56 \cdot 10^9$ |

Table 3: Peak incident rate of the 5 meV detector tubes in McStas and Geant4 simulations, with pyrolytic graphite sample.

3.1.2. Yttrium oxide sample

To give an impression on how the rates change with a different single crystal that is not the worst-case sample but also has a strong Bragg peak, the same simulation and analysis process is repeated using an Y_2O_3 sample.

Fig. 13 demonstrates the change of the time-averaged neutron energy spectrum along the neutrons path at the sample in the scattering characterisation system. The time-averaged neutron intensities acquired by the integration of the energy spectra are shown in Tab. 4.

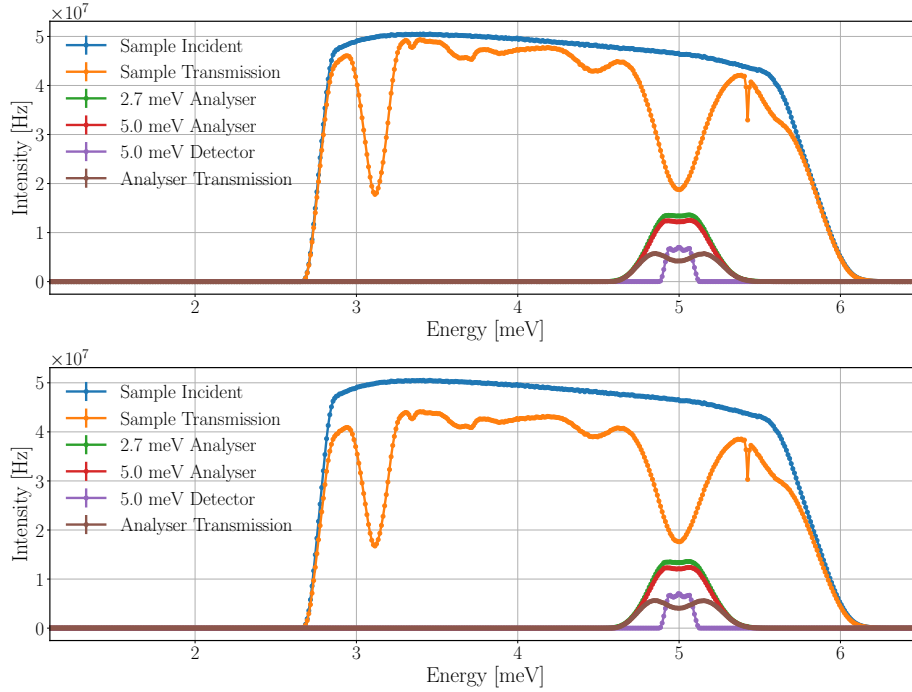


Figure 13: Time-averaged neutron energy spectra at the sample and the scattering characterisation system with an Y_2O_3 single crystal sample in McStas (top) and Geant4 (bottom) simulation. Incident beam on sample (in blue), beam transmitted through the sample (in orange), beam on the set of analysers for 2.7 meV neutrons (in green), beam on the set of analysers for 5.0 meV neutrons (in red), neutrons hitting the detector triplet for 5.0 meV (in purple), beam transmitted through all sets of analysers (in brown). The lines are only joining the points.

The results of the McStas and Geant4 simulations agree with only minor

differences, with the only exception being the transmission spectrum of the sample. The transmission of the sample is 10% higher in McStas, with the same lower intensity bands apparent in the spectrum caused by several crystal planes where the Bragg-criterion is fulfilled for different neutron energies – including the $hkl=-2-2-2$ plane for 5.0 meV. The source of the discrepancy is the absorption process in the sample, as different absorption cross-sections are used in the two simulation tools, as described earlier in section 2.4.

| Position | McStas [Hz] | Geant4 [Hz] |
|-----------------------|----------------------|----------------------|
| Sample Incident | $1.45 \cdot 10^{10}$ | $1.45 \cdot 10^{10}$ |
| Sample Transmission | $1.24 \cdot 10^{10}$ | $1.12 \cdot 10^{10}$ |
| 2.7 meV Analyser | $6.01 \cdot 10^8$ | $6.00 \cdot 10^8$ |
| 5.0 meV Analyser | $5.52 \cdot 10^8$ | $5.48 \cdot 10^8$ |
| 5.0 meV Detector | $1.21 \cdot 10^8$ | $1.21 \cdot 10^8$ |
| Analyser Transmission | $2.96 \cdot 10^8$ | $2.91 \cdot 10^8$ |

Table 4: Time-averaged neutron intensities at the sample and different parts of the scattering characterisation system with an Y_2O_3 single crystal sample in McStas and Geant4 simulation.

The spectrum of the 2.7 meV analyser arc shows that despite the presence of multiple Bragg peaks, it is only the 5.0 meV neutrons that are scattered toward the Q-channel on the $hkl=-2-2-2$ plane. There are, however, three other planes on which the 5.0 meV neutrons are Bragg-scattered ($hkl=2-3-3$, $hkl=1-2-3$ and $hkl=1-3-2$), apparently not toward the Q-channel, causing the slight dip of the peak at 5 meV. The spread of the neutrons after the sample is less significant than it was with pyrolytic graphite, so fewer neutrons are lost on the way toward the 5.0 meV analysers.

The time-averaged neutron intensities in the plane of the 5 meV detector tubes depicted in Fig. 14 show a similar diagonal shape as for pyrolytic graphite because of the energy dependence of the scattering angle both for sample and analysers. The time-averaged incident neutron rate for the detector tubes in Tab. 5 – acquired by integration over the area of the tubes – shows that the

neutrons are distributed more evenly among the tubes as a result of the flattened top of the energy spectrum. The maximum of the time-averaged incident rate for a single tube is found to be 41 MHz.

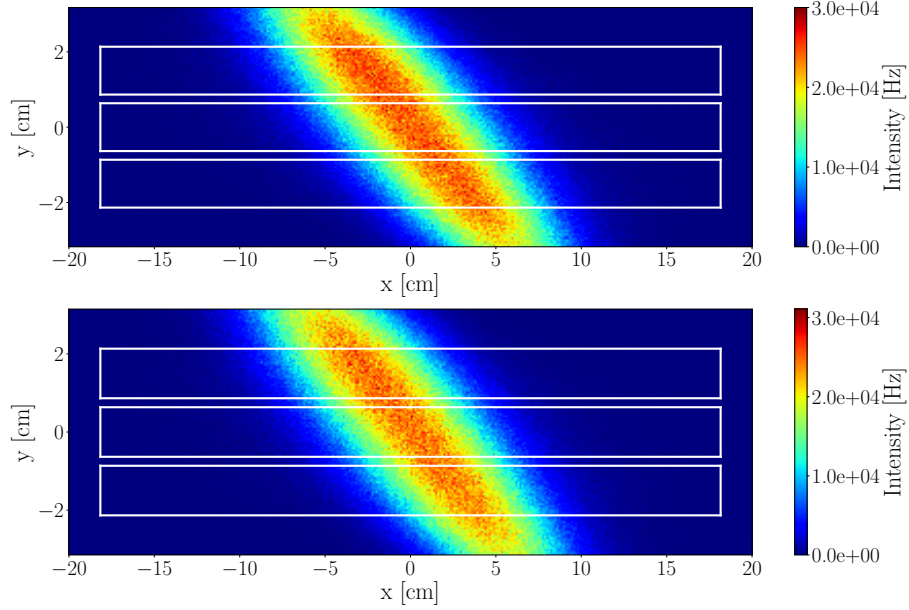


Figure 14: Time-averaged neutron intensities in the plane of the 5 meV detector tubes in McStas (top) and Geant4 (bottom) simulations, with an Y_2O_3 sample. The white lines indicate the outlines of the detector tubes.

Doing the calculations only for the incident neutrons that arrive at the peak of the ToF spectra of each tube, the peak incident rates presented in Tab. 6 are acquired. The peak incident rate for a single detector tube with an Y_2O_3 sample is found to be 1 GHz.

Comparing the time-averaged and peak rates to those acquired for pyrolytic graphite, both are lower by a factor of 1.7 but still on the order of 10 MHz for time-averaged and GHz for peak rates. This means that even with a non-worst-case sample, the rates can be well above the capabilities of the standard ^3He detector tubes.

As demonstrated in this subsection, the McStas and Geant4 simulation results are in excellent agreement, regarding the detector rates. For this reason,

| Detector tube | McStas [Hz] | Geant4 [Hz] |
|---------------|-------------------|-------------------|
| Top | $4.08 \cdot 10^7$ | $4.07 \cdot 10^7$ |
| Central | $4.13 \cdot 10^7$ | $4.13 \cdot 10^7$ |
| Bottom | $3.88 \cdot 10^7$ | $3.86 \cdot 10^7$ |

Table 5: Time-averaged incident neutron rates of the 5 meV detector tubes in McStas and Geant4 simulations, with an Y_2O_3 sample.

| Detector tube | McStas [Hz] | Geant4 [Hz] |
|---------------|-------------------|-------------------|
| Top | $1.03 \cdot 10^9$ | $1.02 \cdot 10^9$ |
| Central | $1.05 \cdot 10^9$ | $1.04 \cdot 10^9$ |
| Bottom | $9.78 \cdot 10^8$ | $9.84 \cdot 10^8$ |

Table 6: Peak incident rates of the 5 meV detector tubes in McStas and Geant4 simulations, with an Y_2O_3 sample.

further simulations are only performed using the Geant4 model of the sample and the scattering characterisation system. In the subsequent subsections multiple parameters are scanned in order to determine their effect on the incident detector rates, and to prove that the results presented above can be regarded as the worst-case incident rates. The change in the incident detector rates due to modifying the studied parameters are expected to have the same trend for all single crystals, so all simulations are done using the Y_2O_3 sample.

3.2. Sample mosaicity

The mosaicity of the sample has multiple effects on the neutron beam Bragg-scattered on a selected scattered plane toward the Q-channel. A sample with higher mosaicity scatters neutrons of a wider energy range, as the higher spread of crystal plane orientations enables them to fulfil the Bragg-criteria. This is also true for a wider incident angle range, meaning that neutrons of a divergent beam with higher incident angle have the possibility to be Bragg-scattered on the selected scattering plane. This higher spread of crystal plane orientations,

on the other hand, lowers the probability of neutrons with energy and incident angle close to the ideal values to be scattered. The cumulative effect is depicted in Fig. 15, showing the energy spectra of the scattered beam at the 2.7 meV analysers and the 5.0 meV detector triplet for an Y_2O_3 sample with different mosaicities.

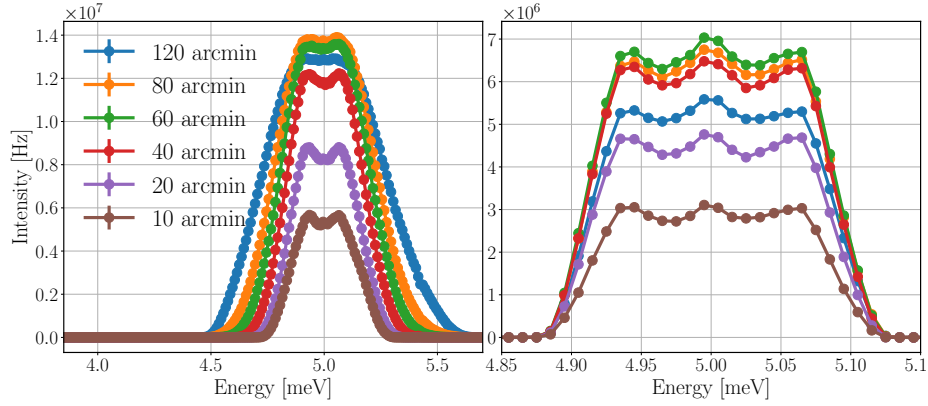


Figure 15: Energy spectra of the neutron beam on the set of analysers for 2.7 meV neutrons (left) and of the neutrons hitting the detector triplet for 5.0 meV (right) for different sample mosaicities with an Y_2O_3 sample. The lines are only joining the points.

The energy spectra of the beam at the 2.7 meV analyser arc shows that the sample with higher mosaicities scatters neutrons of a wider energy range toward the analysers, as expected. It also shows that the intensities at this point of the instrument are getting higher for mosaicities up until 80 arcmin and therefore intensities for 80 arcmin are higher than those for 60 arcmin. It is the energy spectrum of neutrons hitting the detector triplet for 5.0 meV, that shows that these additional neutrons do not reach the detectors, as the highest intensities are found at 60 arcmin sample mosaicity. The resulting time-averaged and peak incident rates of the central detector tube are presented in Tab. 7.

The results are in compliance with the expectation that the highest rates occur when the mosaicity of the sample matches that of the analysers, but also show that within the ± 20 arcmin range it is a less than a 10% effect.

| Mosaicity [arcmin] | Time-averaged rate [Hz] | Peak rate [Hz] |
|--------------------|-------------------------|-------------------|
| 120 | $3.30 \cdot 10^7$ | $8.37 \cdot 10^8$ |
| 80 | $3.99 \cdot 10^7$ | $1.02 \cdot 10^9$ |
| 60 | $4.13 \cdot 10^7$ | $1.04 \cdot 10^9$ |
| 40 | $3.84 \cdot 10^7$ | $9.63 \cdot 10^8$ |
| 20 | $2.80 \cdot 10^7$ | $7.07 \cdot 10^8$ |
| 10 | $1.81 \cdot 10^7$ | $4.54 \cdot 10^8$ |

Table 7: Time-averaged and peak incident neutron rates of the central 5 meV detector tube for different sample mosaicities with an Y_2O_3 sample.

3.3. Analyser mosaicity

The mosaicity of the analysers is a fixed value of 60 arcmin for BIFROST, but it is worth briefly investigating how it would affect the rate of the detector tubes. Fig. 16 depicts the neutron energy spectra of the 5.0 meV detector triplet for different analyser mosaicities with an Y_2O_3 sample with a mosaicity of 60 arcmin.

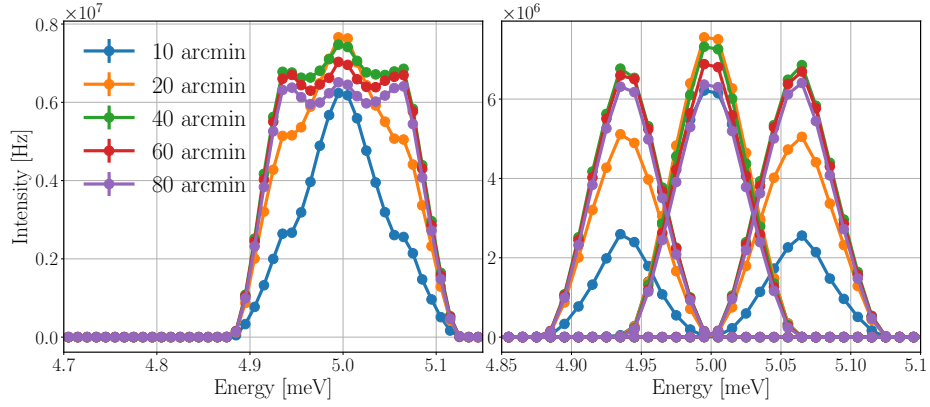


Figure 16: Energy spectra of the neutrons hitting the detector triplet for the three tubes together (left) and separately (right) for 5.0 meV for different analyser mosaicities with an Y_2O_3 sample. The lines are only joining the points.

The spectra of the three tubes separately show that for mosaicities below

| Mosaicity [arcmin] | Time-averaged rate [Hz] | Peak rate [Hz] |
|--------------------|-------------------------|-------------------|
| 10 | $4.10 \cdot 10^7$ | $1.04 \cdot 10^9$ |
| 20 | $4.68 \cdot 10^7$ | $1.19 \cdot 10^9$ |
| 40 | $4.45 \cdot 10^7$ | $1.12 \cdot 10^9$ |
| 60 | $4.13 \cdot 10^7$ | $1.04 \cdot 10^9$ |
| 80 | $3.84 \cdot 10^7$ | $9.64 \cdot 10^8$ |

Table 8: Time-averaged and peak incident neutron rates of the central 5 meV detector tubes for different analyser mosaicities with an Y_2O_3 sample.

40 arcmin the two tubes on the sides are under-illuminated compared to the tube in the centre. In order to apply the prismatic analyser concept, the analyser mosaicity has to be large enough to sufficiently cover all used detectors. Increasing the mosaicity above 40 arcmin, the intensity in all three tubes is getting slightly lower. The resulting time-averaged and peak rates of the central detector tube are presented in Tab. 8.

The results show that the incident rate in a single detector tube could be 13–14% higher in the central tube with 20 arcmin mosaicity compared to the result with 60 arcmin, but the mosaicity has to be higher to apply the prismatic analyser concept, and in the range of 40–80 arcmin the change is less than 10%.

3.4. Sample size

Sample size is the limiting factor in many scientific cases, as it is not easy to grow large samples of some types. The beam delivery system of BIFROST is optimised for sample cross-sections up to $15 \times 15 \text{ mm}^2$ but the realistic sample sizes for the intended applications are much smaller than that, with an expected minimum sample size going down to 1 mm^3 . In order to investigate the effect of the sample size on the incident detector rates, cylindrical samples with equal diameter and height are simulated. The energy spectra of the scattered beam at the 2.7 meV analysers and the 5.0 meV detector triplet for an Y_2O_3 sample of different sizes are depicted in Fig. 17. The resulting time-averaged and peak

incident rates of the central detector tube are presented in Tab. 9.

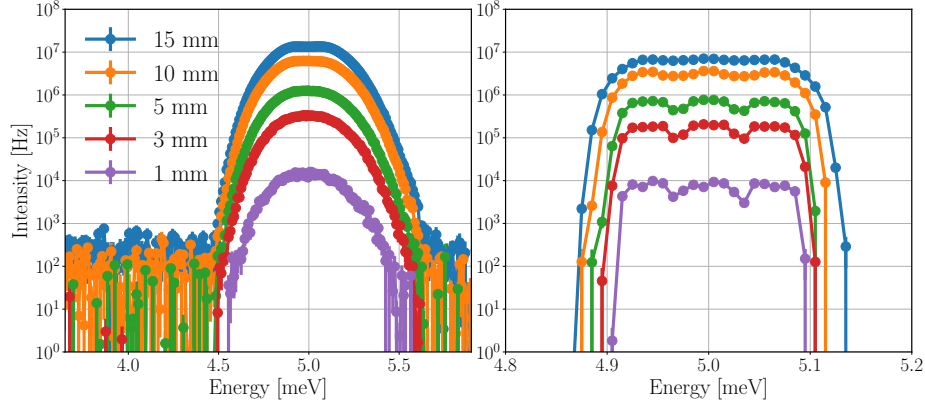


Figure 17: Energy spectra of the neutron beam on the set of analysers for 2.7 meV neutrons (left) and of the neutrons hitting the detector triplet for 5.0 meV (right) for an Y_2O_3 sample of different sizes. The diameter and height of the cylindrical samples are equal, with the magnitude indicated in the legend. The lines are only joining the points.

| Sample size [mm] | Time-averaged rate [Hz] | Peak rate [Hz] |
|------------------|-------------------------|-------------------|
| 15 | $4.13 \cdot 10^7$ | $1.04 \cdot 10^9$ |
| 10 | $1.93 \cdot 10^7$ | $4.93 \cdot 10^8$ |
| 5 | $3.92 \cdot 10^6$ | $1.01 \cdot 10^8$ |
| 3 | $1.04 \cdot 10^6$ | $2.8 \cdot 10^7$ |
| 1 | $4.6 \cdot 10^4$ | $1.6 \cdot 10^6$ |

Table 9: Time-averaged and peak incident neutron rates of the central 5 meV detector tube for an Y_2O_3 sample of different sizes.

As expected, the larger the sample, the higher the intensities are. By reducing the sample size parameter (height and diameter) from 15 mm to 5 mm and 1 mm, the time-averaged incident rate of the center tube drops by a factor of 10.5 and 900 respectively. Due to the better resolutions in case of smaller samples the drop in the peak incident rate is lower, a factor of 10.3 for 5 mm and a factor of 650 for 1 mm.

Another effect of the better resolution is visible in the energy spectra of the detector triplets, where the three-peak structure is more apparent for smaller samples.

3.5. Pulse-shaping chopper opening time

The energy resolution of the instrument can be increased at the cost of neutron intensity by modifying the opening time of the pulse-shaping chopper. The flux on the sample and the detectors are both expected to drop significantly in case of the high resolution setting of BIFROST, where the pulse-shaping chopper opening time is merely 0.1 ms, compared to the high flux mode achieved by an opening time of 5 ms.

The energy spectra of the neutron beam at the sample and the 5.0 meV detector triplet for an Y_2O_3 sample for different pulse-shaping chopper opening times are depicted in Fig. 18, with the time-averaged intensities presented in Tab. 10. The resulting time-averaged and peak rates of the central detector tube are presented in Tab. 11.

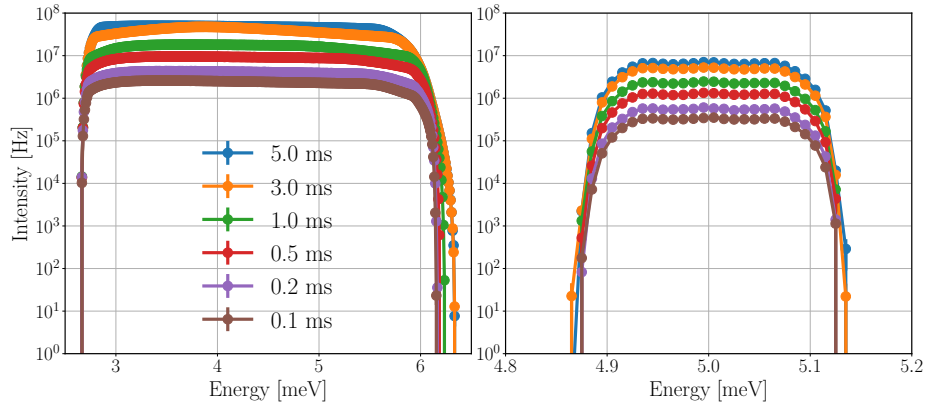


Figure 18: Energy spectra of the neutron beam on the sample (left) and of the neutrons hitting the detector triplet for 5.0 meV (right) for an Y_2O_3 sample for different pulse-shaping chopper opening times. The lines are only joining the points.

As expected, the time-averaged rates on the sample and on the detectors decrease with shorter pulse-shaping chopper opening times. The difference in

| PSC opening time [ms] | Sample [Hz] | 5.0 meV detectors [Hz] |
|-----------------------|----------------------|------------------------|
| 5.0 | $1.45 \cdot 10^{10}$ | $1.21 \cdot 10^8$ |
| 3.0 | $1.19 \cdot 10^{10}$ | $9.08 \cdot 10^7$ |
| 1.0 | $5.09 \cdot 10^9$ | $4.21 \cdot 10^7$ |
| 0.5 | $2.76 \cdot 10^9$ | $2.28 \cdot 10^7$ |
| 0.2 | $1.26 \cdot 10^9$ | $1.02 \cdot 10^7$ |
| 0.1 | $7.47 \cdot 10^8$ | $5.99 \cdot 10^6$ |

Table 10: Time-averaged neutron intensities at the sample and the 5.0 meV detector tubes with an Y_2O_3 sample for different pulse-shaping chopper (PSC) opening times.

| PSC opening time [ms] | Time-averaged rate [Hz] | Peak rate [Hz] |
|-----------------------|-------------------------|-------------------|
| 5.0 | $4.13 \cdot 10^7$ | $1.04 \cdot 10^9$ |
| 3.0 | $3.11 \cdot 10^7$ | $1.04 \cdot 10^9$ |
| 1.0 | $1.45 \cdot 10^7$ | $8.41 \cdot 10^8$ |
| 0.5 | $7.83 \cdot 10^6$ | $5.33 \cdot 10^8$ |
| 0.2 | $3.50 \cdot 10^6$ | $2.57 \cdot 10^8$ |
| 0.1 | $2.05 \cdot 10^6$ | $1.53 \cdot 10^8$ |

Table 11: Time-averaged and peak incident neutron rates of the central 5 meV detector tube with Y_2O_3 sample for different pulse-shaping chopper opening times.

time-averaged incident rates between the high flux mode (5 ms) and the high resolution mode (0.1 ms) is approximately a factor of 20 for both the detector triplet and the central tube. Regarding the peak rates in the central tube, however, this drop is less apparent. The highest rate for a 3 ms opening time is the same (within statistical uncertainty) as the rate for 5 ms, and the difference compared to 0.1 ms opening time is only a factor of 6.8. The reason for this difference is the better ToF resolution with shorter pulse-shaping chopper opening times. The higher time-averaged rates are distributed over a longer period of time on the detectors, as it is demonstrated in Fig. 19, showing the effect of the pulse-shaping chopper opening time on the ToF spectrum of neutrons hitting the detectors.

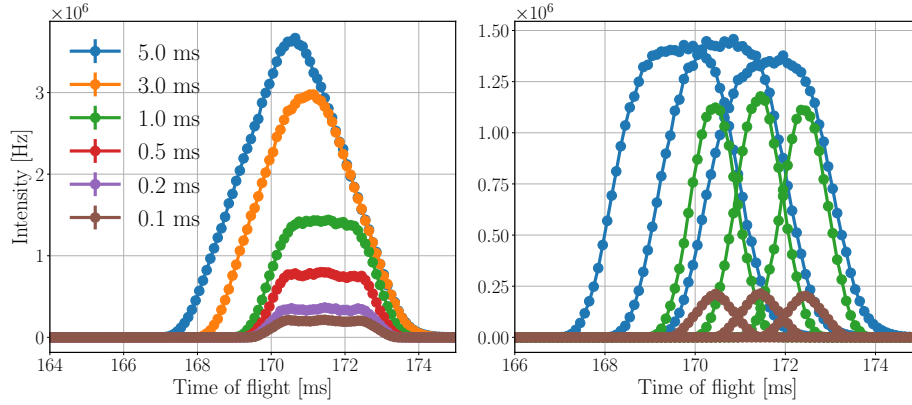


Figure 19: ToF spectrum of neutrons at the 5 meV detector triplet for the three tubes together (on left) and separately (on right) with Y_2O_3 sample for different pulse-shaping chopper opening times. The separated spectrum is not displayed in all cases in the righthand figure to avoid the figure being overcrowded. The lines are only joining the points.

The longer the opening time, the broader the ToF peaks are. This directly affects the energy resolution, and also increases the dead time in case of saturation. If one tube of a triplet is saturated, then none of the three can read out data, as they are connected in series. This means that for the detector triplet in the presented case for 5.0 ms pulse-shaping chopper opening time no data is recorded for more than 7 ms.

4. Elastic peak rates in representative operational conditions

Parameters chosen in section 3.1 correspond to possible worst-case scenarios and rates acquired are far above the capabilities of ^3He tubes. However, the combination of a strongly scattering large sample and the highest flux mode is rather artificial, so it is worth evaluating a more representative operational scenario.

BIFROST is designed for small samples, as sample size is the limiting factor in many science cases. Hence, cm-size crystals are not to be expected very often, only large samples with small magnetic moments, and therefore small magnetic Bragg peak intensity. There is another parameter directly affecting the intensities but not discussed yet, the accelerator power of the ESS source. As mentioned earlier in section 2.2, the source power of 5 MW is used for the simulations. That is the eventual operational power of ESS, however, it will initially operate at 2 MW. The intensities are expected to scale linearly with the source power, and in terms of the simulations, the source power is in fact only a multiplication factor in the ESS butterfly moderator source component.

For these reasons, the following parameters are selected to define the rates in a more representative operational case: 2 MW source power, 1 ms pulse-shaping chopper opening time, an Y_2O_3 single crystal sample with a height and diameter of 3 mm, and mosaicity of 60 arcmin. The time-averaged energy spectra of the neutron beam at the sample and different positions of the scattering characterisation system, acquired with these parameters are demonstrated in Fig. 20, along with the integral values presented in Tab. 12. The resulting time-averaged and peak incident neutron rates of the 5 meV detector tubes are presented in Tab. 13.

The combined effect of the lower source power, shorter pulse-shaping chopper opening time and smaller sample size (cross-section) decreased the time-averaged neutron intensity on the sample significantly, by a factor of 136 compared to the worst-case scenario with an Y_2O_3 sample. Due to the reduced sample thickness, the transmission through the sample is increased to 94% from

| Position | Intensity [Hz] |
|-----------------------|-------------------|
| Sample Incident | $1.07 \cdot 10^8$ |
| Sample Transmission | $1.01 \cdot 10^8$ |
| 2.7 meV Analyser | $1.84 \cdot 10^6$ |
| 5.0 meV Analyser | $1.70 \cdot 10^6$ |
| 5.0 meV Detector | $4.1 \cdot 10^5$ |
| Analyser Transmission | $8.71 \cdot 10^5$ |

Table 12: Time-averaged neutron intensities at the sample and different parts of the scattering characterisation system with an Y_2O_3 sample, using 2 MW source power, 1 ms PSC opening time, 3 mm height and diameter sample size and 60 arcmin sample mosaicity.

| Detector tube | Time-averaged rate [Hz] | Peak rate [Hz] |
|---------------|-------------------------|------------------|
| Top | $1.3 \cdot 10^5$ | $9.3 \cdot 10^6$ |
| Central | $1.5 \cdot 10^5$ | $9.9 \cdot 10^6$ |
| Bottom | $1.3 \cdot 10^5$ | $9.4 \cdot 10^6$ |

Table 13: Time-averaged and peak incident neutron rates of the 5 meV detector tubes with an Y_2O_3 sample, using 2 MW source power, 1 ms PSC opening time, 3 mm height and diameter sample size and 60 arcmin sample mosaicity.

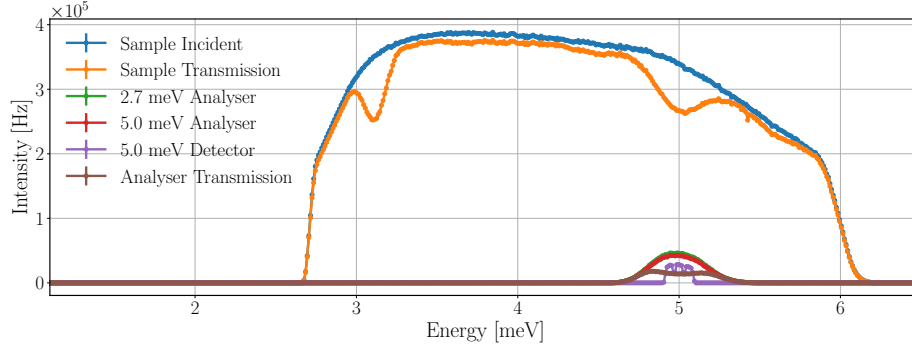


Figure 20: Time-averaged neutron energy spectra at the sample and the scattering characterisation system with an Y_2O_3 sample in Geant4 simulation using 2 MW source power, 3 mm height and diameter sample size and 60 arcmin sample mosaicity. Incident beam on sample (in blue), beam transmitted through the sample (in orange), beam on the set of analysers for 2.7 meV neutrons (in green), beam on the set of analysers for 5.0 meV neutrons (in red), neutrons hitting the detector triplet for 5.0 meV (in purple), beam transmitted through all sets of analysers (in brown). The lines are only joining the points.

77%, as a result of lower absorption and weaker Bragg peaks. The lower incident intensity on the sample, and weaker Bragg peak lead to a drop by a factor of 322–326 in the time-averaged neutron intensity on both the 2.7 and 5.0 meV analysers. The drop in the time-averaged neutron intensity on the 5.0 meV detector triplet is a little bit lower, a factor of 295, due to the smaller divergence and better energy resolution of the neutron beam compared to the worst-case scenario. The highest time-averaged incident neutron rate on a single tube is found to be 0.15 MHz, that means a drop by a factor of 275, but the peak incident rate is 9.9 MHz, which is lower only by a factor 105 compared to the worst-case results. The numbers and reduction factors are in accordance with previous simulations in sections 3.4 and 3.5 where the effect of sample size and pulse-shaping chopper opening time were investigated separately.

There are factors not taken into account in the current study, that may further reduce the intensities slightly on the detectors, like the non-ideal transmission of the filtering system, and the effect of the divergence jaws applicable for reducing the angular spread of neutrons.

5. Simulation with full scattering characterisation system

The previous sections were aimed to define the highest incident rates a detector tube can experience using different instrument and sample parameters in case of a coherent elastic (Bragg) peak. For this purpose different single crystal samples were used with the simulation model including only a single Q-channel, that was rotated according to the Bragg-angle to get the highest possible rate for the detectors. This section intends to present the incident detector rates in case of incoherent elastic peaks with a standard calibration sample, and to demonstrate the use of the full simulation model of the scattering characterisation system with all nine Q-channels.

The sample selected for this simulation is vanadium, that is assumed to be an incoherent elastic scatterer which scatters isotropically and therefore it is used to calibrate the incident neutron intensity and the detector efficiencies in neutron spectrometers [34]. As mentioned earlier, in previous simulations the single Q-channel present in the model was rotated according to the Bragg-angle of the sample for 5 meV neutrons ($\Theta = 37.067^\circ$ for pyrolytic graphite and $\Theta = 41.169^\circ$ for Y_2O_3). For vanadium the rotation of the nine Q-channels is arbitrarily selected in a way to have $2\Theta = 90^\circ$ scattering angle for the central Q-channel, as depicted in Fig. 21.

The instrument and sample parameters are the same as for the worst-case scenario: 5 MW source power, 5 ms pulse-shaping chopper opening time, 15 mm sample height and diameter. The time-averaged energy spectra of the neutron beam at the sample and different positions of the central Q-channel are demonstrated in Fig. 22, along with the integral values presented in Tab. 14.

The transmission through the sample is only 23% – much lower than it is for pyrolytic graphite (80%) or Y_2O_3 (77%) – with no peaks missing from the spectrum, as expected from a sample scattering mainly incoherently. The wide spectrum of the 2.7 meV analysers also shows that neutrons are not coming from an elastic peak, but despite the wider energy range, a logarithmic scale is needed as the integrated intensity is more than 2 orders of magnitude lower

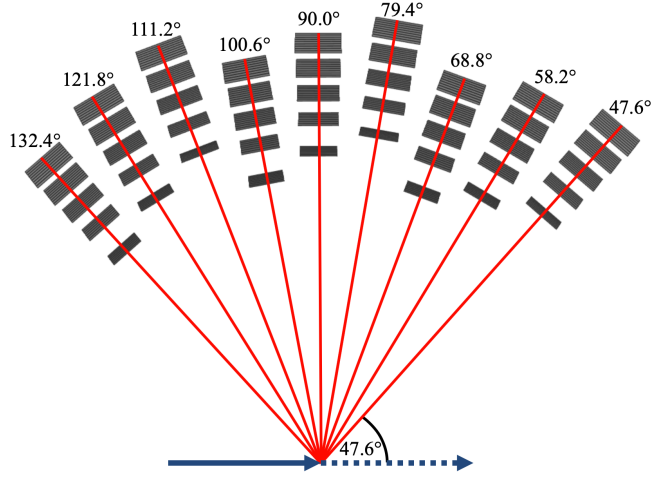


Figure 21: Top view schematic figure of the scattering characterisation system model with all nine Q-channels. The red lines and corresponding angles indicate the scattering angle for the centre of each Q-channel.

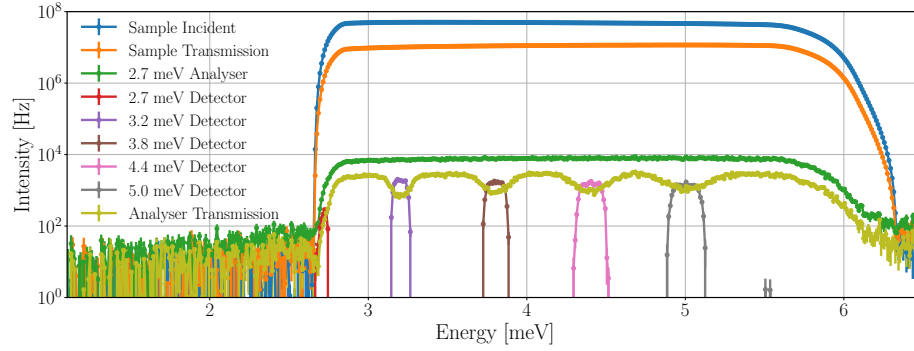


Figure 22: Time-averaged neutron energy spectra at the sample and in the central Q-channel with a vanadium sample in Geant4 simulation. Incident beam on sample (in blue), beam transmitted through the sample (in orange), beam on the set of analysers for 2.7 meV neutrons (in green), neutrons hitting the detector triplets for energies 2.7–5.0 meV (in red–grey), beam transmitted through all sets of analysers (in mustard). The lines are only joining the points.

than experienced with previous samples. The spectrum of the beam transmitted through all five sets of analysers clearly shows neutrons missing because they are selected by the analysers. These neutrons appear in the spectra of the detectors, that show that peaks are narrower for lower energies as a result of the better

| Position | Intensity [Hz] |
|-----------------------|--------------------------|
| Sample Incident | $1.45 \cdot 10^{10}$ |
| Sample Transmission | $3.36 \cdot 10^9$ |
| 2.7 meV Analyser | $2.55 \cdot 10^6$ |
| 2.7 meV Detector | $1.6 \pm 0.3 \cdot 10^3$ |
| 3.2 meV Detector | $1.7 \pm 0.1 \cdot 10^4$ |
| 3.8 meV Detector | $2.1 \pm 0.1 \cdot 10^4$ |
| 4.4 meV Detector | $2.4 \pm 0.1 \cdot 10^4$ |
| 5.0 meV Detector | $2.6 \pm 0.2 \cdot 10^4$ |
| Analyser Transmission | $7.20 \cdot 10^5$ |

Table 14: Time-averaged neutron intensities at the sample and different parts of the scattering characterisation system with a vanadium sample in Geant4 simulation.

energy resolution of the analysers for lower energies. The incident intensity on the 2.7 meV detectors is much lower than in other detectors due the energy range selected by the bandwidth chopper. The resulting time-averaged incident rates are higher for higher energies, with the maximum of 26 ± 2 kHz for the 5.0 meV detector triplet. For a single detector tube the highest rates are found for the central 5.0 meV detector with a time-averaged intensity of 9 ± 1 kHz and peak intensity of 0.3 ± 0.1 MHz. The time-averaged incident neutron rates of all detector triplets in all Q-channels are presented in Tab. 15.

The trends in the results demonstrate the combination of three effects. In each Q-channel the detector triplets for higher energies experience higher incident rates due to the wider energy ranges selected by the analysers, as shown for the central Q-channel earlier in this section. The second effect has roots in the “triple stagger” geometry and the asymmetry of the Q-channels described in section 2.1. The sample–analyser distances in Q-channels 1, 4, 7 are shorter, and in channels 3, 6, 9 are longer than distances in channels 2, 5, 8. The shorter distances increase the rates visibly because neutrons are not collimated by Bragg-scattering on the sample and therefore their spread at longer dis-

| Q-channel | Scattering angle [°] | 2.7 meV [kHz] | 3.2 meV [kHz] | 3.8 meV [kHz] | 4.4 meV [kHz] | 5.0 meV [kHz] |
|-----------|----------------------|---------------|---------------|---------------|---------------|---------------|
| 1 | 132.4 | 2.8±0.3 | 23±1 | 28±1 | 30±2 | 31±2 |
| 2 | 121.8 | 2.2±0.3 | 20±1 | 23±1 | 27±2 | 30±2 |
| 3 | 111.2 | 1.9±0.3 | 18±1 | 21±1 | 24±1 | 25±1 |
| 4 | 100.6 | 2.3±0.3 | 21±1 | 25±1 | 28±2 | 30±2 |
| 5 | 90.0 | 1.6±0.3 | 17±1 | 21±1 | 24±1 | 26±2 |
| 6 | 79.4 | 1.4±0.2 | 15±1 | 18±1 | 21±1 | 22±2 |
| 7 | 68.8 | 2.0±0.2 | 18±1 | 22±1 | 25±2 | 26±1 |
| 8 | 58.2 | 1.4±0.2 | 15±1 | 19±1 | 22±1 | 23±1 |
| 9 | 47.6 | 1.1±0.2 | 13±1 | 16±1 | 19±1 | 21±1 |

Table 15: Time-averaged neutron intensities of the 5 detector triplets in all nine Q-channels with a vanadium sample.

tances becomes important. This effect on the rates is somewhat blurred by the third effect caused by the anisotropy of the scattering cross-section in vanadium. For the observed energies the scattering cross-section of vanadium is higher for higher scattering angles [35]. This results in generally higher rates for Q-channels positioned for higher scattering angles, but due to the asymmetry of the adjacent Q-channels, it is most apparent when comparing Q-channels of the same symmetry, like 2, 5, 8.

6. Conclusions

BIFROST is an indirect ToF spectrometer at ESS and one of the first eight instruments to be constructed. One of the most challenging aspects of its operation is the rate capability and in particular the peak instantaneous rate capability, i.e. the number of neutrons hitting the detector per channel at the peak of the neutron pulse. There is no scientific reason to measure the elastic peaks as they are considered background for this instrument, however it is vital that the detectors are not degraded by such intensity and remain capable of measuring weak inelastic signals, as soon as possible after saturation. This implies that the detector aspects of recovery time and high rate tolerance have to be carefully evaluated by measurements to prove that scientific performance will be intact.

A detailed methodology for acquiring the results is presented. The full simulation of the instrument from source to detector position is carried out using multiple simulation software packages. A flexible model of the sample and the scattering characterisation system of BIFROST is implemented in both McStas and Geant4 and a comparison of their strengths and weaknesses is presented. The capability of both simulation tools is enhanced by the NCrystal library and associated tools. The first application of the special NCrystal pyrolytic graphite is presented, demonstrating its capabilities for modelling analysers for neutron scattering applications.

McStas is capable of simulating instruments as long as 160 m, and even handling beam splitting to some extent, to treat simulations with multiple set of analysers, however the latter comes with great complexity and some limitations, as it is not within the natural usage of this simulation software. Geant4 on the other hand is not suited to simulate the beam transport system of an instrument, but with the use of NCrystal, it is an entirely appropriate tool for a scattering characterisation system with any level of geometrical complexity, and even offers the possibility to include parts like filtering system and cross-talk shielding, and taking into account back-scattering. The results of the McStas and Geant4 model of the scattering characterisation system are compared us-

ing various single crystal samples. The results show perfect agreement with the only exception being the transmission through the sample where a difference less than 10% is found in one case, due to the more detailed modelling of absorption in Geant4.

With this knowledge at hand a choice was made to combine the McStas model of the beam transport system and the Geant4 model of the sample and the scattering characterisation system using the MCPL tools. Using this model, the impact of sample type, sample and analyser mosaicity, sample size, and pulse shaping chopper opening time is studied on the incident detector rates.

The primary purpose of this work is to provide an estimate of the incident rates that are anticipated for the BIFROST instrument, and also to demonstrate the use of powerful simulation tools for the correct treatment of neutron transport in crystalline materials. The incident detector rates anticipated for the BIFROST instrument for different configurations are presented. For instrument configurations and sample parameters representing worst-case conditions, it is determined that the peak rate can reach the value of 1–1.7 GHz for a single detector tube with time-averaged rates of 40–70 MHz. These tubes are expected to reach saturation well below that, at 50–100 kHz. These tubes will also be saturated for a minimum of 5 ms, but the saturation deadtime for detecting signals is more like 7 ms because the counting detector tubes are coupled in triplets.

To overcome challenges caused by these rates, an operational evaluation of a measurement strategy will be the key to the successful operation of this instrument. More “everyday” realistic samples give a lower rate challenge, however these samples will still saturate detectors.

A simulation with the full analyser system is presented using a common calibration sample. This model can now be used to predict experimental conditions from specific proposed samples, i.e. sample size and composition for experiment planning purposes for users.

The results here show the potential power of source to detector simulation for neutron scattering. These simulations are now possible due to tools recently

developed. It is now possible to realistically simulate very complex systems.

Acknowledgements

This work has been supported by the In-Kind collaboration between ESS ERIC (contract number: NIK5.4 #10 [ESS]) and the Hungarian Academy of Sciences, Centre for Energy Research (MTA EK). The authors would like to thank the DMSC Computing Centre which provided the computing resources for the simulations [36]. Richard Hall-Wilton, Kalliopi Kanaki and Thomas Kittelmann would like to acknowledge support from BrightnESS [EU Horizon2020 grant 676548].

References

- [1] S. Peggs, et al., ESS Technical Design Report, ESS 2013-001, <http://esss.se/scientific-technical-documentation> (2013).
- [2] R. Garoby, et al., The european spallation source design, Physica Scripta 93 (1) (2018) 014001.
URL <http://stacks.iop.org/1402-4896/93/i=1/a=014001>
- [3] K. Andersen, et al., The instrument suite of the european spallation source, Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment 957 (2020) 163402. doi:<https://doi.org/10.1016/j.nima.2020.163402>.
URL <http://www.sciencedirect.com/science/article/pii/S0168900220300097>
- [4] H. Ronnow, K. Lefmann, N. B. Christensen, C. Niedermayer, F. Juranyi, M. Marko, J. O. Birk, M. Bertelsen, Bifrost instrument proposal.
URL https://ess-public-legacy.esss.se/sites/default/files/bifrost_proposal_may_2014.pdf

- [5] P. G. Freeman, et al., CAMEA ESS The continuous angle multi-energy analysis indirect geometry spectrometer for the European Spallation Source, EPJ Web Conf. 83 (2015) 03005. [arXiv:1406.5945](#), doi: 10.1051/epjconf/20158303005.
- [6] G. F. Knoll, Radiation detection and measurement; 4th ed., Wiley, New York, NY, 2010.
URL <https://cds.cern.ch/record/1300754>
- [7] K. Kanaki, et al., Detector rates for the Small Angle Neutron Scattering instruments at the European Spallation Source, Journal of Instrumentation 13 (07) (2018) P07016.
URL <http://stacks.iop.org/1748-0221/13/i=07/a=P07016>
- [8] I. Stefanescu, M. Christensen, J. Fenske, R. Hall-Wilton, P. Henry, O. Kirstein, M. Mller, G. Nowak, D. Pooley, D. Raspino, N. Rhodes, J. Šaroun, J. Schefer, E. Schooneveld, J. Sykora, W. Schweika, Neutron detectors for the ESS diffractometers, Journal of Instrumentation 12 (01) (2017) P01019–P01019. doi:10.1088/1748-0221/12/01/p01019.
URL <https://doi.org/10.1088%2F1748-0221%2F12%2F01%2Fp01019>
- [9] F. Piscitelli and others, The Multi-Blade Boron-10-based neutron detector for high intensity neutron reflectometry at ESS, JINST 12 (2017) 03013.
- [10] K. Kanaki, et al., Simulation tools for detector and instrument design, Physica B: Condensed Matter 551 (2018) 386 – 389, the 11th International Conference on Neutron Scattering (ICNS 2017). doi:<https://doi.org/10.1016/j.physb.2018.03.025>.
URL <http://www.sciencedirect.com/science/article/pii/S0921452618302175>
- [11] I. Stefanescu, M. Christensen, R. Hall-Wilton, S. Holm-Dahlin, K. Iversen, M. Klein, D. Mannix, J. Schefer, C. Schmidt, W. Schweika, U. Stuhr, Performance study of the jalousie detector baseline design for the ESS

- thermal powder diffractometer HEIMDAL through GEANT4 simulations, *Journal of Instrumentation* 14 (10) (2019) P10020–P10020. doi:10.1088/1748-0221/14/10/p10020.
URL <https://doi.org/10.1088%2F1748-0221%2F14%2F10%2Fp10020>
- [12] X.-X. Cai, T. Kittelmann, Ncrystal: A library for thermal neutron transport, *Computer Physics Communications* 246 (2020) 106851. doi:<https://doi.org/10.1016/j.cpc.2019.07.015>.
URL <http://www.sciencedirect.com/science/article/pii/S0010465519302280>
- [13] T. Kittelmann and others, Monte Carlo Particle Lists: MCPL, *Computer Physics Communications* 218 (2017) 17–42. doi:10.1016/j.cpc.2017.04.012.
- [14] <https://mctools.github.io/mcpl/>, MCPL documentation and GitHub repository (2017).
- [15] K. Lefmann, K. Nielsen, Mcstas, a general software package for neutron ray-tracing simulations, *Neutron News* 10 (3) (1999) 20–23. doi:10.1080/10448639908233684.
- [16] P. Willendrup, E. Farhi, K. Lefmann, McStas 1.7 - a new version of the flexible Monte Carlo neutron scattering package, *Physica B: Condensed Matter* 350 (13, Supplement) (2004) E735 – E737, proceedings of the Third European Conference on Neutron Scattering. doi:10.1016/j.physb.2004.03.193.
- [17] S. Agostinelli, et al., GEANT4: A Simulation toolkit, *Nucl. Instrum. Meth. A* 506 (2003) 250–303. doi:10.1016/S0168-9002(03)01368-8.
- [18] J. Allison, et al., Geant4 developments and applications, *IEEE Trans. Nucl. Sci.* 53 (2006) 270. doi:10.1109/TNS.2006.869826.
- [19] J. Allison, et al., Recent developments in Geant4, *Nucl. Instrum. Meth. A* doi:10.1016/j.nima.2016.06.125.

- [20] M. Marko, F. Groitl, J. Birk, P. Freeman, K. Lefmann, N. Christensen, C. Niedermayer, F. Juranyi, J. Lass, A. Hansen, H. Ronnow, Prototype of the novel camea concept – a backend for neutron spectrometers, *Review of Scientific Instruments* 89 (2018) 015105. doi:10.1063/1.5018233.
- [21] F. Groitl, D. Graf, J. O. Birk, M. Mark, M. Bartkowiak, U. Filges, C. Niedermayer, C. Regg, H. M. Rnnow, Camea – a novel multiplexing analyzer for neutron spectroscopy, *Review of Scientific Instruments* 87 (3) (2016) 035109. arXiv:<https://aip.scitation.org/doi/pdf/10.1063/1.4943208>, doi:10.1063/1.4943208.
URL <https://aip.scitation.org/doi/abs/10.1063/1.4943208>
- [22] S. Holm-Dahlin, M. A. Olsen, M. Bertelsen, J. O. Birk, K. Lefmann, Optimization of performance, price, and background of long neutron guides for european spallation source, *Quantum Beam Science* 3 (3) (2019) 16. doi:10.3390/qubs3030016.
URL <http://dx.doi.org/10.3390/qubs3030016>
- [23] European Spallation Source ESS BIFROST Instrument, <https://europeanspallationsource.se/instruments/bifrost>.
- [24] F. Groitl, E. Rantsiou, M. Bartkowiak, U. Filges, D. Graf, C. Niedermayer, C. Regg, H. M. Rnnow, A combined radial collimator and cooled beryllium filter for neutron scattering, *Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment* 819 (2016) 99 – 103. doi:<https://doi.org/10.1016/j.nima.2016.02.056>.
URL <http://www.sciencedirect.com/science/article/pii/S0168900216002321>
- [25] D. F. R. Mildner, M. Arif, S. A. Werner, Neutron transmission through pyrolytic graphite monochromators, *Journal of Applied Crystallography* 34 (3) (2001) 258–262. doi:10.1107/S0021889801002254.
URL <https://doi.org/10.1107/S0021889801002254>

- [26] J. O. Birk, M. Mark, P. G. Freeman, J. Jacobsen, R. L. Hansen, N. B. Christensen, C. Niedermayer, M. Mnsson, H. M. Rnnow, K. Lefmann, Prismatic analyser concept for neutron spectrometers, Review of Scientific Instruments 85 (11) (2014) 113908. [arXiv:https://doi.org/10.1063/1.4901160](https://doi.org/10.1063/1.4901160), doi:10.1063/1.4901160.
URL <https://doi.org/10.1063/1.4901160>
- [27] E. Knudsen, A. Prodi, J. Baltser, M. Thomsen, P. Willendrup, M. Snchez del Ro, C. Ferrero, E. Farhi, K. Haldrup, A. Vickery, R. Feidenhans'l, K. Mortensen, M. Nielsen, H. Poulsen, S. Schmidt, K. Lefmann, Mcxtrace: A monte carlo software package for simulating x-ray optics, beamlines and experiments, Journal of Applied Crystallography 46 (2013) 679–696. doi:10.1107/S0021889813007991.
- [28] C.J. Werner (editor), MCNP Users Manual - Code Version 6.2, Tech. Rep. LA-UR-17-29981, Los Alamos National Laboratory (2017).
- [29] BIFROST McStas model.
URL <https://bitbucket.org/europeanspallationsource/nosg-baselines/src/dev/BIFROST/Optics/McStas/>
- [30] K. H. Andersen, M. Bertelsen, L. Zanini, E. B. Klinkby, T. Schönfeldt, P. M. Bentley, J. Saroun, Optimization of moderators and beam extraction at the ESS, Journal of Applied Crystallography 51 (2) (2018) 264–281. doi:10.1107/S1600576718002406.
URL <https://doi.org/10.1107/S1600576718002406>
- [31] L. Zanini, K. Batkov, E. Klinkby, F. Mezei, T. Schnfeldt, A. Takibayev, The neutron moderators for the european spallation source, Journal of Physics: Conference Series 1021 (2018) 012066. doi:10.1088/1742-6596/1021/1/012066.
URL <https://doi.org/10.1088%2F1742-6596%2F1021%2F1%2F012066>
- [32] P. Willendrup, E. Farhi, E. Knudsen, U. Figles, K. Lefmann, J. Stein, User

and programmers guide to the neutron ray-tracing package mcstas, version 2.6 (2020) P72–P73.

- [33] T. Kittelmann, et al., Geant4 based simulations for novel neutron detector development, *J. Phys: Conf. Ser.* 513 (2014) 022017. doi:10.1088/1742-6596/513/2/022017.
- [34] J. Mayers, The use of vanadium as a scattering standard for pulsed source neutron spectrometers, *Nuclear Instruments and Methods in Physics Research* 221 (3) (1984) 609 – 618. doi:[https://doi.org/10.1016/0167-5087\(84\)90073-5](https://doi.org/10.1016/0167-5087(84)90073-5).
URL <http://www.sciencedirect.com/science/article/pii/0167508784900735>
- [35] J. Mayers, Contribution of inelastic scattering to the vanadium differential scattering cross section; implications for the calibration of neutron spectrometers, *Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment* 281 (3) (1989) 654 – 656. doi:[https://doi.org/10.1016/0168-9002\(89\)91504-0](https://doi.org/10.1016/0168-9002(89)91504-0).
URL <http://www.sciencedirect.com/science/article/pii/0168900289915040>
- [36] Data Management & Software Centre webpage.
URL <https://europeanspallationsource.se/data-management-software/computing-centre>