# Imaging and Identification of Single Nanoplastic Particles and Agglomerates

Ambika Shorny<sup>1</sup>, Fritz Steiner<sup>1</sup>, Helmut Hörner<sup>1</sup>, and Sarah M. Skoff<sup>1,\*</sup>


<sup>1</sup>Atominstitut, Technische Universität Wien, Stadionallee 2, Vienna, 1020, Austria \*sarah.skoff@tuwien.ac.at

## ABSTRACT


Pollution by nanoplastic is a growing environmental and health concern. Currently the extent of nanoplastic in the environment can only be cumbersomely and indirectly estimated but not measured. To be able to quantify the extent of the problem, detection methods that can identify nanoplastic particles that are smaller than 1  $\mu$ m are critically needed. Here, we employ surface-enhanced Raman scattering (SERS) to image and identify single nanoplastic particles down to 100 nm in size. We can differentiate between single particles and agglomerates and our method allows an improvement in detection speed of  $10^7$  compared to state-of-the art surface-enhanced Raman imaging. Being able to resolve single particles allows to measure the SERS enhancement factor on individual nanoplastic particles instead of averaging over a concentration without spatial information. Our results thus contribute to the better understanding and employment of SERS for nanoplastic detection and present an important step for the development of future sensors.

The global increase in plastic production and disposal has resulted in large amounts of plastic that end up in our environment. This was even further intensified by the recent COVID19 pandemic which required the widespread increased use of single-use-plastics for testing and personal protection equipment<sup>1</sup>. Fragmentation of this plastic waste in the environment leads to microand nanoplastics which are dispersed even more easily and are harder to contain. Nanoplastics, commonly defined as particles with one dimension being smaller than 1  $\mu$ m<sup>2</sup>, are in particular very difficult to detect because of their small size but can cause great harm as these particles are able to cross the blood-brain barrier in living organisms and even enter the placenta<sup>3</sup>. Harmful effects on human health that arise due to continued exposure to micro- and especially nanoplastic are therefore an active field of research<sup>4,5</sup>.

As a first step in tackling the nanoplastic problem, suitable detection methods are required that can visualize and quantify the problem at hand. Thus far, the larger microplastic is mostly detected by visual inspection, Fourier transform infrared (FTIR) spectroscopy<sup>6</sup> and Raman microscopy<sup>7</sup>. While visual inspection<sup>8</sup> is fastest and most widely used, it depends heavily on the interpretation of the observer and also lacks chemical information. Further, it cannot be applied for nanoparticles. FTIR and Raman spectroscopy are able to provide information on the vibrational structure of the detected molecules but still have some disadvantages. FTIR is restricted to larger particle sizes on the order of 20  $\mu$ m and suffers from strong background signals from water which requires samples to be carefully dried. Raman microscopy works for smaller particles and can be

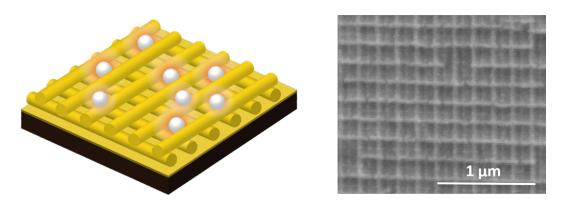

**Figure 1.** Sketch of SERS substrate with polystyrene beads (left) and scanning electron microscope image of SERS substrate used for the detection of nanoplastic (right).

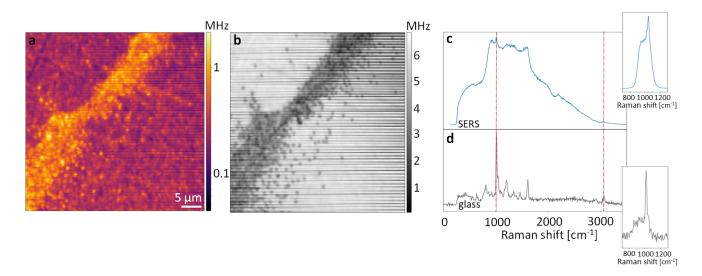

**Figure 2.** Agglomerates of and single 300 nm PS particles on the SERS substrate visualised (a) via surface-enhanced Raman imaging and (b) via recording the backscattered laser beam. (c) Raman spectra obtained from a single nanoparticle on a SERS substrate and on a (d) reference glass substrate. Two unbroadened Raman modes that are visible in both spectra are indicated via dashed lines. The insets shows the respective Raman spectra with the bandpass filter to select the brightest C-C stretch vibrational mode.

used for spectroscopy in water but Raman scattering is an inherently weak process, with only about one in 10<sup>8</sup> photons of the excitation beam contributing to the Raman signal<sup>9</sup>. Other methods such as thermal analysis yields chemical information but is destructive and information about number and size is not provided<sup>10</sup> and scanning electron microscopy (SEM) can provide information on the geometry of small particles but requires careful sample preparation and makes it difficult to obtain chemical information. For these reasons, methods are required that offer good signal levels to ensure fast detection times, while still yielding information on the type of the analyte.

Here, we show that by employing surface-enhanced Raman scattering (SERS), the Raman signals can be vastly enhanced by more than three orders of magnitude and by combining SERS with confocal microscopy, we are able to image and identify single plastic particles down to sizes of 100 nm. While surface-enhanced Raman spectroscopy has shown great benefits in the detection of small molecular coverages, where higher enhancement factors can be obtained as molecules can sample the smallest gaps of nanostructures<sup>11–13</sup>, it is only very recently that it has been applied for the detection of nanoparticles and aerosols<sup>14–17</sup> and in particular nanoplastics<sup>18–27</sup>. With our system, we are able to image and identify single nanoplastic particles and obtain the SERS enhancement factor without the need to average over different concentrations. Additionally, we demonstrate that we can decrease the image acquisition time for spatially-resolved nanoplastic detection by filtering the strongest Raman band instead of using Raman mapping and simultaneously imaging the particles by recording the scattered laser light, which we demonstrate to result in an improvement by a factor of 10<sup>7</sup> for spatially-resolved nanoplastic detection. Our detection method is not limited to any particular nanoplastic concentration and thus provides an important step towards the development of nanoplastic sensors in the future.

#### 1 Materials and Methods

#### 1.1 SERS substrate

SERS enhancement is generally implemented in two different ways: either the analyte is mixed with a solution of plasmonic nanoparticles and then deposited on a non-SERS substrate<sup>20, 22, 24, 26, 28</sup> or a plasmonic substrate is manufactured beforehand and the analyte is then deposited on it<sup>18, 21, 23, 25, 27</sup>. The latter case has the advantage that in principle a higher reproduceability of the results should be possible. We therefore follow the latter method by employing a nano-patterned gold substrate (Pico Foundry<sup>29</sup>) which consists of 3D crossed gold nanowires that form a square grid with a side length of 100 nm (Fig. 1). In our case, the enhancement of the Raman signal is due to electromagnetic enhancement. This means that as gold has a plasmon resonance in the visible, localized surface plasmons get excited and thus light gets tightly confined by the gold nanostructure forming "hot spots". The scattering of particles that sample these "hot spots" then also gets intensified by the plasmonic enhancement.

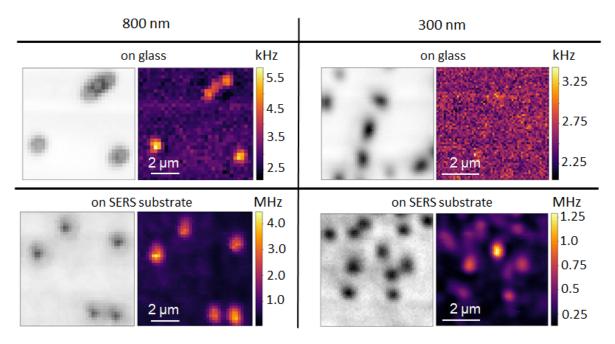


Figure 3. Comparison of confocal images taken of 800 nm and 300 nm polystyrene particles on a glass and a SERS substrate.

#### **1.2 Experimental Procedures**


The nanoplastic used was purchased from Sigma Aldrich in three different sizes: 800 nm, 300 nm and 100 nm diameter spherical polystyrene (PS) particles. They come dispersed in water (10%) and before depositing them on the different substrates we dilute them further (1:1300) with ultrapure water to be able to also obtain single nanoparticles on the substrates. For detection of the nanoplastic on either substrate, a homebuilt confocal microscope is used. A 500 mW laser at 638 nm is connected to the microscope setup via a fiber and the excitation beam is focussed via a Mitutoyo 50x objective with an NA of 0.55 onto the sample, which can be scanned by a stepper motor and a piezo in closed- or open-loop operation. The maximum power at the sample is about 24 mW. The signal is collected by the same objective and transmitted through a dichroic mirror and 75  $\mu$ m pinhole. It is split into two beams by a 50/50 beamsplitter, which are filtered and then coupled into two optical fibers which lead to two single photon counting modules (SPCM). One of the beams is attenuated by neutral density filters before entering the detector, while for the other beam, we filter out the excitation beam completely by two longpass filters (LP 650 nm) and pick out the brightest part of the Raman band of polystyrene by a bandpass filter (BP 680 nm).

The backscattered light field consists of a contribution which is reflected by the substrate and a contribution that is affected by any particle present on the substrate by being either scattered or absorbed. When this light field is imaged as a function of position via the confocal microscope, the particles thus appear as shadows (Fig. 2 (b)), whereas in the other arm, where the brightest Raman mode is selected those particles that are plastic are bright on a dark background (Fig. 2 (a)). The maximum scan range of the stepper motor than can move the sample is 25 x 25 mm and that of the piezo is  $100 \times 100 \,\mu$ m. Typical pixel sizes of a piezo scan are between 100-300 nm. The typical time taken for closed-loop operation is 0.2 s per pixel, where 0.1 s integration time is used and the rest of the time is taken by the program to stabilize the piezo at the required position. To get an overview of a big area, the pixel size can be increased and the piezo can be operated without feedback with a lower integration time of down to 10 ms/pixel. Concerning our observed Raman scattering rate of the particles on the SERS substrate, integration times of less than 1 ms/pixel are possible and thus at the moment we are only limited by the way our piezo stage is implemented in our experimental control.

To obtain a wide Raman spectrum from about 290-4000  $\text{cm}^{-1}$ , the bandpass filter is removed and the Raman emission can be sent to a spectrometer (Shamrock SR-303i, Andor Technology). We have used two different gratings of 300 l/mm and 600 l/mm and typical integration times are between 1-10 s.

## 2 Results and Discussion


An example of an image of 300 nm polystyrene beads is displayed in Fig. 2 (a) and (b), where polystyrene beads with a diameter of 300 nm are imaged on the SERS substrate over an area of  $40 \times 40 \,\mu$ m. Sending the Raman scattered light at the position of the particles onto a spectrometer yields the known spectrum of polystyrene. Fig. 2 shows the obtained spectrum of

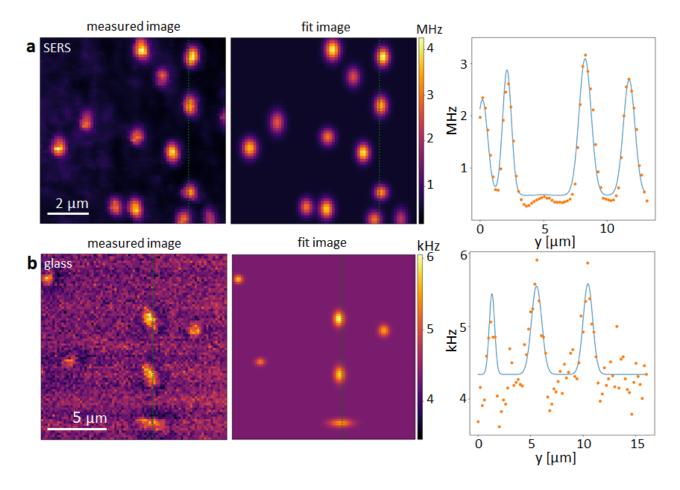

**Figure 4.** Measured Raman images and their corresponding fitted images for (a) the SERS substrate and (b) a glass substrate for 800 nm PS particles. On the right of these images a 2D slice at the position indicated in the respective images is plotted to further visualise the gain in absolute signal and in signal to background ratio when imaging these nanoparticles via SERS.

an individual particle on a (c) SERS substrate and (d) a reference glass substrate. When the particles are placed on the SERS target some Raman modes get broadened<sup>30</sup> (Fig. 2 (c)) but one C-C stretch mode ( $v_{24}$ ) and one C-H stretch mode ( $v_7$ ) remain unaffected and thus can still be used for identifying the particles as polystyrene<sup>31</sup>.

To obtain fast Raman images, we forego Raman mapping by filtering the strongest Raman mode at 1003 cm<sup>-1</sup> (C-C stretch mode) by a bandpass filter and sending this light to an SPCM to create the images. The spectrum with the bandpass filter introduced is displayed in the insets of Fig. 2 (c) and (d). The surface-enhanced Raman image displayed in Fig. 2 (a) shows agglomerates and individual 300 nm PS particles. Here, a log scale is applied to the colorbar to best visualise single particles as well as the brighter agglomerate. This stresses the importance of also being able to image individual nanoparticles to make sure the detection method is really suitable for the desired size range and down to the level of single particle detection. To compare the results obtained from the SERS substrate to the Raman images obtained on a glass substrate, Fig. 3 shows confocal images taken of 800 nm and 300 nm polystyrene particles on both types of substrates. For 800 nm particles, both conventional Raman and SERS can resolve the particles, albeit SERS at a better signal to noise ratio. Importantly, however, for the 300 nm particles, the conventional Raman signal is already very small and thus the particles are barely visible, whereas on the SERS substrate they can be clearly distinguished from the background (Fig. 3). This clearly shows that our here introduced SERS platform enables Raman detection of small nanoplastic particles in a facile fashion. By directly using appropriate filters instead of employing Raman mapping, these surface-enhanced Raman images can also be obtained at a much greater speed and integration times of < 1ms/pixel are feasible (Supplementary Information and Figure S2). To evaluate the experimental enhancement factor of the SERS substrate and quantify the improvement of the images in terms of contrast, we fit the images to the point spread function of the scatterers (Fig. 4), which is well approximated by a two-dimensional Gaussian<sup>32</sup>. This yields

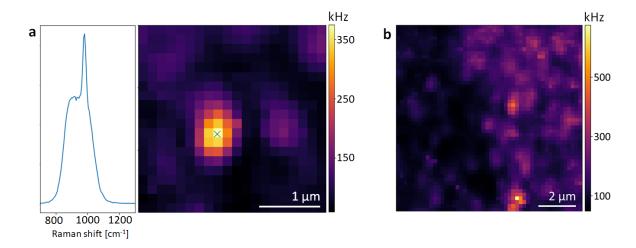

**Figure 5.** (a) Spectrum and surface-enhanced Raman image of a 100 nm PS particle. (b) Raman image of agglomerates and single 100 nm nanoplastic particles

the peak intensity of the Raman emission and we calculate the experimental enhancement factor as:

$$EF = \frac{I_{SERS}}{I_{glass}},$$
(1)

where I<sub>SERS</sub> is the peak intensity of the particles on the SERS substrate and I<sub>glass</sub> is the peak intensity on the glass substrate for the same Raman band. In both cases the background has been substracted. We find that the SERS substrate enhances the Raman signal by more than 3 orders of magnitude (Table 1). Such a high enhancement factor has not been measured yet on single nanoplastic particles<sup>25</sup> and our results show that being able to image single particles is key to obtaining the true enhancement factor for detection of particles of such small sizes. Employing SERS, we also increase the signal to background ratio by more than an order of magnitude. This is visualised in Fig. 4, where a slice through an image is taken, depicting the particles as peaks in this 2D plot. The values for the enhancement factor and the signal to background ratio are given in table 1 for a laser power of  $23.32 \pm 0.85$  mW at the target, which was the highest power used and the best in terms of contrast for the glass substrate. While we are only able to see Raman images of the plastic particles on glass for excitation laser powers of > 4 mW, for the SERS substrate powers as low as 150  $\mu$ W show a clear signal. For lower powers the experimentally evaluated enhancement factor therefore represents a lower limit. In terms of signal to background ratio, there is a tendency that even for the SERS substrate this ratio decreases as the size of the particles decreases. This is not surprising as the size of the particles gets smaller and thus the fraction of light from our diffraction limited excitation laser beam of about 1.45  $\mu$ m that illuminates the particles with respect to the illuminated area of the background, decreases as well. Employing SERS for nanoparticles is generally a very new field that has only recently shown promising results<sup>18, 20–28</sup>, but previously, enhancement factors for SERS have been determined for different concentrations of the analyte without any spatial information, which leaves a large uncertainty on the number of particles that actually contribute to the signal<sup>20</sup>. Often, the advantage of SERS is illustrated by stating a minimum detectable concentration of nanoplastic<sup>22,23,26</sup>, which is still larger than the concentrations of hundreds to thousands of ng/l that are expected in the environment<sup>33–35</sup> and also lacks spatial information. Being able to image single particles is thus an important step towards developing precise sensors in the future that are able to detect low concentrations of single nanoplastic particles rather than only agglomerates and for understanding and optimizing SERS microscopy for nanoparticles in general. Due to the enhancement in Raman signal obtained from the SERS substrate we are able to detect and image individual particles with a diameter of only 100 nm, which is far below the diffraction limited spot of our laser beam. A SERS image of a 100

particle size	S/B glass	S/B SERS substrate	EF
800	$1.849 \pm 0.098$	$6.425\pm0.839$	979 ±189
300	$1.116\pm0.033$	$3.743 \pm 1.165$	$2219 \pm 1109$
100	-	$2.205 {\pm}~0.327$	-

**Table 1.** Experimental enhancement factor (EF) and the signal to background ratio (S/B) for single nanoparticles of different sizes on a glass or SERS substrate.

nm particle with the corresponding spectrum with bandpass filter can be seen in Fig. 5(a) and a larger image of agglomerates of and single 100 nm particles is seen in Fig. 5(b). This clearly shows that the plasmonic enhancement allows for imaging nanoplastic particles that otherwise remain invisible.

## 3 Conclusion

In a time where not much is known yet about the distribution of nanoplastic particles in the environment<sup>36</sup>, techniques need to be developed to quantify and monitor this pollution. This is crucial as it is expected that nanoplastic is at least as widespread as microplastics, while having the potential to be more harmful to living organisms<sup>37, 38</sup>.

In this study, we have demonstrated the imaging and identification of single nanoplastic particles and of agglomerates down to sizes of 100 nm via surface-enhanced Raman microscopy with a gain in acquisition speed of 10<sup>7</sup> compared to state-of-the art micro-Raman measurements<sup>20, 22, 23, 26</sup>. We also demonstrate the variation in signal strength between single particles and agglomerates and hence the need for spatial information to assess the detection limits of the developed methods. For particles of 300 nm and 800 nm, we compare our measurements to Raman microscopy on a non-SERS substrate and find enhancement factors of more than three orders of magnitude. To our knowledge, these are the largest enhancement factors seen so far for individually imaged nanoplastic particles<sup>25</sup> and are needed to image particles of such a small size in such a facile way.

Instead of using Raman mapping we employ filtering of the brightest Raman band of the analyte to obain direct Raman images using our confocal microscope and simultaneously image the particles by recording the backscattered laser light. Compared to Raman mapping which usually takes 1-10 s integration time per pixel, this is a vast improvement<sup>18</sup> as integration times of less than 1 ms/pixel become feasible for an image of individual nanoplastic particles. Our way of imaging can be even further improved by employing more narrowband filters<sup>43</sup> possibly also for different Raman bands. In addition, if one wanted to record the whole spectrum for a wide wavelength range, spectra only need to be taken for those coordinates where bright particles are visible in the confocal images. Currently, simultaneous imaging of the backscattered laser light lets us image particles down to sizes of 300 nm, but this resolution can be further improved by better control over the phase of the laser beam<sup>44</sup>. While not being particle specific, this way of imaging also allows a quick overview of interesting parts of the sample with minimal laser powers.

Nanoplastic concentrations in the environment span a wide range<sup>33–35</sup> and hence it is crucial to develop sensitive methods that also work on the single particle level. Our results not only aid in understanding and employing SERS enhancement for nanoplastic detection but are an important step for the process of developing the first in-situ detector to monitor nanoplastics in the environment in the future.

### 4 Data availability

The datasets used and/or analysed during the current study are available from the corresponding author on reasonable request.

#### **5** Acknowledgments

We want to thank A. Rauschenbeutel for the loan of equipment. This work was supported by the Austrian Research Promotion Agency (FFG) (no. 884447, PhoQus2D). S.M.S. acknowledges funding from the Austrian Science Fund (FWF) (no. V934-N, Quantoom).

## 6 Authors' contributions

A.S. performed all the measurements and analysed the data. F.S. analysed the data and H.H. set up the experimental control. S.M.S. conceived the idea, supervised the project and wrote the manuscript with input from all the authors.

#### 7 Competing interests

The authors declare no competing interests.

#### References

- 1. Adyel, T. M. Accumulation of plastic waste during COVID-19. *Science* **369**, 1314–1315, DOI: 10.1126/science.abd9925 (2020). Publisher: American Association for the Advancement of Science.
- 2. Cai, H. *et al.* Analysis of environmental nanoplastics: Progress and challenges. *Chem. Eng. J.* 410, 128208, DOI: 10.1016/j.cej.2020.128208 (2021).

- 3. Campanale, Massarelli, Savino, Locaputo & Uricchio. A Detailed Review Study on Potential Effects of Microplastics and Additives of Concern on Human Health. *Int. J. Environ. Res. Public Heal.* 17, 1212, DOI: 10.3390/ijerph17041212 (2020).
- 4. Gruber, E. S. *et al.* To Waste or Not to Waste: Questioning Potential Health Risks of Micro- and Nanoplastics with a Focus on Their Ingestion and Potential Carcinogenicity. *Expo. Heal.* DOI: 10.1007/s12403-022-00470-8 (2022).
- Liu, X. *et al.* Bioeffects of Inhaled Nanoplastics on Neurons and Alteration of Animal Behaviors through Deposition in the Brain. *Nano Lett.* 22, 1091–1099, DOI: 10.1021/acs.nanolett.1c04184 (2022). Publisher: American Chemical Society.
- 6. Käppler, A. *et al.* Analysis of environmental microplastics by vibrational microspectroscopy: FTIR, Raman or both? *Anal. Bioanal. Chem.* **408**, 8377–8391, DOI: 10.1007/s00216-016-9956-3 (2016).
- Hale, R. C., Seeley, M. E., La Guardia, M. J., Mai, L. & Zeng, E. Y. A Global Perspective on Microplastics. J. Geophys. Res. Ocean. 125, DOI: 10.1029/2018JC014719 (2020).
- 8. Ivleva, N. P., Wiesheu, A. C. & Niessner, R. Microplastic in Aquatic Ecosystems. *Angewandte Chemie Int. Ed.* 56, 1720–1739, DOI: 10.1002/anie.201606957 (2017).
- 9. Gardiner, D. J. & Graves, P. R. (eds.) *Practical Raman Spectroscopy* (Springer Berlin Heidelberg, Berlin, Heidelberg, 1989).
- Shim, W. J., Hong, S. H. & Eo, S. E. Identification methods in microplastic analysis: a review. *Anal. Methods* 9, 1384–1391, DOI: 10.1039/C6AY02558G (2017).
- Kim, J. *et al.* Synthesis and Single-Particle Surface-Enhanced Raman Scattering Study of Plasmonic Tripod Nanoframes with Y-Shaped Hot-Zones. *Nano Lett.* 20, 4362–4369, DOI: 10.1021/acs.nanolett.0c01084 (2020). Publisher: American Chemical Society.
- Demirtaş, O., Doğanay, D., Murat Öztürk, I., Emrah Unalan, H. & Bek, A. Facile preparation of nanoparticle based SERS substrates for trace molecule detection. *Phys. Chem. Chem. Phys.* 22, 21139–21146, DOI: 10.1039/D0CP01866J (2020). Publisher: Royal Society of Chemistry.
- Zhang, L. *et al.* Large-scale growth of sharp gold nano-cones for single-molecule SERS detection. *RSC Adv.* 6, 2882–2887, DOI: 10.1039/C5RA22321K (2016). Publisher: The Royal Society of Chemistry.
- Fu, Y. *et al.* Surface-Enhanced Raman Spectroscopy: A Facile and Rapid Method for the Chemical Component Study of Individual Atmospheric Aerosol. *Environ. Sci. & Technol.* 51, 6260–6267, DOI: 10.1021/acs.est.6b05910 (2017). Publisher: American Chemical Society.
- **15.** Sivaprakasam, V. & Hart, M. B. Surface-Enhanced Raman Spectroscopy for Environmental Monitoring of Aerosols. *ACS Omega* **6**, 10150–10159, DOI: 10.1021/acsomega.1c00207 (2021). Publisher: American Chemical Society.
- Yoo, H., Lee, H., Park, C., Shin, D. & Ro, C.-U. Novel Single-Particle Analytical Technique for Submicron Atmospheric Aerosols: Combined Use of Dark-Field Scattering and Surface-Enhanced Raman Spectroscopy. *Anal. Chem.* 94, 13028– 13035, DOI: 10.1021/acs.analchem.2c01696 (2022). Publisher: American Chemical Society.
- Ali Tahir, M. *et al.* Klarite as a label-free SERS-based assay: a promising approach for atmospheric bioaerosol detection. *Analyst* 145, 277–285, DOI: 10.1039/C9AN01715A (2020). Publisher: Royal Society of Chemistry.
- 18. Chang, L. *et al.* Nanowell-enhanced Raman spectroscopy enables the visualization and quantification of nanoplastics in the environment. *Environ. Sci. Nano* 9, 542–553, DOI: 10.1039/D1EN00945A (2022). Publisher: Royal Society of Chemistry.
- **19.** Fang, C. *et al.* Identification and visualisation of microplastics/ nanoplastics by Raman imaging (ii): Smaller than the diffraction limit of laser? *Water Res.* **183**, 116046, DOI: 10.1016/j.watres.2020.116046 (2020).
- Hu, R., Zhang, K., Wang, W., Wei, L. & Lai, Y. Quantitative and sensitive analysis of polystyrene nanoplastics down to 50 nm by surface-enhanced Raman spectroscopy in water. *J. Hazard. Mater.* 429, 128388, DOI: 10.1016/j.jhazmat.2022. 128388 (2022).
- 21. Jeon, Y. *et al.* Detection of nanoplastics based on surface-enhanced Raman scattering with silver nanowire arrays on regenerated cellulose films. *Carbohydr. Polym.* 272, 118470, DOI: 10.1016/j.carbpol.2021.118470 (2021).
- Kihara, S. *et al.* Detecting polystyrene nanoplastics using filter paper-based surface-enhanced Raman spectroscopy. *RSC Adv.* 12, 20519–20522, DOI: 10.1039/D2RA03395J (2022). Publisher: Royal Society of Chemistry.
- Lê, Q. T. *et al.* Nanostructured Raman substrates for the sensitive detection of submicrometer-sized plastic pollutants in water. *J. Hazard. Mater.* 402, 123499, DOI: 10.1016/j.jhazmat.2020.123499 (2021).
- 24. Lv, L. *et al.* In situ surface-enhanced Raman spectroscopy for detecting microplastics and nanoplastics in aquatic environments. *Sci. The Total. Environ.* 728, 138449, DOI: 10.1016/j.scitotenv.2020.138449 (2020).

- 25. Xu, G. *et al.* Surface-Enhanced Raman Spectroscopy Facilitates the Detection of Microplastics <1 μm in the Environment. *Environ. Sci. & Technol.* 54, 15594–15603, DOI: 10.1021/acs.est.0c02317 (2020). Publisher: American Chemical Society.
- 26. Zhou, X.-X., Liu, R., Hao, L.-T. & Liu, J.-F. Identification of polystyrene nanoplastics using surface enhanced Raman spectroscopy. *Talanta* 221, 121552, DOI: 10.1016/j.talanta.2020.121552 (2021).
- 27. Yang, Q. *et al.* Identification of Trace Polystyrene Nanoplastics Down to 50 nm by the Hyphenated Method of Filtration and Surface-Enhanced Raman Spectroscopy Based on Silver Nanowire Membranes. *Environ. Sci. & Technol.* DOI: 10.1021/acs.est.2c02584 (2022). Publisher: American Chemical Society.
- 28. Lin, P.-Y., Wu, I.-H., Tsai, C.-Y., Kirankumar, R. & Hsieh, S. Detecting the release of plastic particles in packaged drinking water under simulated light irradiation using surface-enhanced Raman spectroscopy. *Anal. Chimica Acta* 1198, 339516, DOI: 10.1016/j.aca.2022.339516 (2022).
- 29. PiCO SERS.
- **30.** Anema, J. R., Brolo, A. G., Felten, A. & Bittencourt, C. Surface-enhanced Raman scattering from polystyrene on gold clusters. *J. Raman Spectrosc.* **41**, 745–751, DOI: 10.1002/jrs.2504 (2010). Publisher: John Wiley & Sons, Ltd.
- **31.** Sears, W. M., Hunt, J. L. & Stevens, J. R. Raman scattering from polymerizing styrene. I. Vibrational mode analysis. *The J. Chem. Phys.* **75**, 1589–1598, DOI: 10.1063/1.442262 (1981). Publisher: American Institute of Physics.
- 32. Qu, X., Wu, D., Mets, L. & Scherer, N. F. Nanometer-localized multiple single-molecule fluorescence microscopy. *Proc. Natl. Acad. Sci.* 101, 11298–11303, DOI: 10.1073/pnas.0402155101 (2004). Publisher: Proceedings of the National Academy of Sciences.
- 33. Li, Z. *et al.* Impacts of polystyrene nanoplastics at the environmentally relevant and sub-lethal concentrations on the oxidative stress, immune responses, and gut microbiota to grass carp (Ctenopharyngodon idella). *J. Hazard. Mater.* 441, 129995, DOI: 10.1016/j.jhazmat.2022.129995 (2023).
- **34.** Enders, K., Lenz, R., Stedmon, C. A. & Nielsen, T. G. Abundance, size and polymer composition of marine microplastics ≥10 µm in the Atlantic Ocean and their modelled vertical distribution. *Mar. Pollut. Bull.* **100**, 70–81, DOI: 10.1016/j. marpolbul.2015.09.027 (2015).
- **35.** Pabortsava, K. & Lampitt, R. S. High concentrations of plastic hidden beneath the surface of the Atlantic Ocean. *Nat. Commun.* **11**, 4073, DOI: 10.1038/s41467-020-17932-9 (2020). Number: 1 Publisher: Nature Publishing Group.
- 36. Lins, T. F., O'Brien, A. M., Zargartalebi, M. & Sinton, D. Nanoplastic State and Fate in Aquatic Environments: Multiscale Modeling. *Environ. Sci. & Technol.* 56, 4017–4028, DOI: 10.1021/acs.est.1c03922 (2022). Publisher: American Chemical Society.
- 37. Lehner, R., Weder, C., Petri-Fink, A. & Rothen-Rutishauser, B. Emergence of Nanoplastic in the Environment and Possible Impact on Human Health. *Environ. Sci. & Technol.* 53, 1748–1765, DOI: 10.1021/acs.est.8b05512 (2019). Publisher: American Chemical Society.
- **38.** Liu, L. *et al.* Cellular internalization and release of polystyrene microplastics and nanoplastics. *Sci. The Total. Environ.* **779**, 146523, DOI: 10.1016/j.scitotenv.2021.146523 (2021).
- Xie, L., Gong, K., Liu, Y. & Zhang, L. Strategies and Challenges of Identifying Nanoplastics in Environment by Surface-Enhanced Raman Spectroscopy. *Environ. Sci. & Technol.* 57, 25–43, DOI: 10.1021/acs.est.2c07416 (2023). Publisher: American Chemical Society.
- 40. Pérez, A., Prada, Y. A., Cabanzo, R., González, C. I. & Mejía-Ospino, E. Diagnosis of chagas disease from human blood serum using surface-enhanced Raman scattering (SERS) spectroscopy and chemometric methods. *Sens. Bio-Sensing Res.* 21, 40–45, DOI: 10.1016/j.sbsr.2018.10.003 (2018).
- **41.** Xu, M. *et al.* Internalization and toxicity: A preliminary study of effects of nanoplastic particles on human lung epithelial cell. *Sci. The Total. Environ.* **694**, 133794, DOI: 10.1016/j.scitotenv.2019.133794 (2019).
- **42.** Goul, R. *et al.* Quantitative analysis of surface enhanced Raman spectroscopy of Rhodamine 6G using a composite graphene and plasmonic Au nanoparticle substrate. *Carbon* **111**, 386–392, DOI: 10.1016/j.carbon.2016.10.019 (2017).
- **43.** Sun, C. *et al.* Tunable narrow-band single-channel add-drop integrated optical filter with ultrawide FSR. *PhotoniX* **3**, 12, DOI: 10.1186/s43074-022-00056-2 (2022).
- 44. Taylor, R. W. & Sandoghdar, V. Interferometric Scattering Microscopy: Seeing Single Nanoparticles and Molecules via Rayleigh Scattering. *Nano Lett.* **19**, 4827–4835, DOI: 10.1021/acs.nanolett.9b01822 (2019).