

Exploring Novel Object Recognition and Spontaneous Location Recognition Analysis Techniques in Alzheimer’s Mice

Soham Bafana¹, Radha Raghuraman^{2, 3}, and S. Abid Hussaini^{2, 3}

²Taub Institute for Research on Alzheimer’s disease and the Aging Brain, Columbia University Irving Medical Center

³Department of Pathology and Cell Biology

December 13, 2023

1 Abstract

Understanding object recognition patterns in mice is crucial for advancing behavioral neuroscience and has significant implications for human health, particularly in the realm of Alzheimer’s research. This study is centered on the development, application, and evaluation of a state-of-the-art computational pipeline designed to analyze such behaviors, specifically focusing on Novel Object Recognition (NOR) and Spontaneous Location Recognition (SLR) tasks. The pipeline integrates three advanced computational models—Any-Maze for initial data collection, DeepLabCut for detailed pose estimation, and Convolutional Neural Networks (CNNs) for nuanced behavioral classification. Employed across four distinct mouse groups, this pipeline demonstrated high levels of accuracy and robustness. Despite certain challenges like video quality limitations and the need for manual calculations, the results affirm the pipeline’s efficacy and potential for scalability. The study serves as a proof of concept for a multidimensional computational approach to behavioral neuroscience, emphasizing the pipeline’s versatility and readiness for future, more complex analyses.

2 Introduction

Understanding intricate cognitive behaviors in animal models, especially mice, has gained paramount importance in behavioral neuroscience. This significance is accentuated by the rise of computational models that offer unprecedented precision and depth in behavioral analysis. Our study focuses on leveraging this computational power to explore Novel Object Recognition (NOR) and Spontaneous Location Recognition (SLR) in mice, both of which serve as critical platforms for examining cognitive function.

NOR and SLR are widely recognized as invaluable tools for probing various aspects of cognitive function, including memory and spatial recognition. They have been employed in the study of neurodegenerative conditions such as Alzheimer’s disease (AD), providing crucial insights into cognitive impairments associated with these disorders [6]. NOR stands out as a robust, non-invasive method for assessing memory and recognition. Its quick execution time and minimal stress on subjects make it particularly attractive for behavioral analysis [1]. The NOR task comprises three phases: habituation, familiarization, and testing, each designed to provide unique insights into cognitive function. Intriguingly, the hippocampus and perirhinal cortex play pivotal roles in NOR, contributing to different aspects of object recognition and offering a fertile ground for exploring how computational models can dissect these complex neural interactions [1, 3].

On the other hand, SLR provides a nuanced platform for measuring spatial pattern separation abilities and has also been utilized to study cognitive impairments in neurodegenerative disorders [7]. What sets SLR apart is its flexibility; researchers can manipulate object configurations in the

arena to challenge pattern separation processes and offer deeper insights into spatial memory and cognition [7].

Central to our investigation are three cutting-edge technologies: Any-Maze, DeepLabCut, and Convolutional Neural Networks (CNNs). Any-Maze enables us to define zones around objects and quantify the time spent exploring them, allowing for nuanced behavioral analysis. DeepLabCut, leveraging deep learning, offers markerless pose estimation to track and quantify mice exploration of objects, providing a dynamic perspective on behavioral patterns. Meanwhile, CNNs, utilizing the EXPLORE[5] open-source tools, allow for complex video analysis to detect and segment objects, thereby enabling a comprehensive understanding of exploration behavior.

The synergy of these computational models aims to provide a multidimensional analysis framework specifically tailored for NOR and SLR. This approach promises to revolutionize the precision and depth of behavioral analyses, setting the stage for future research directions in understanding cognitive impairments, particularly in neurodegenerative conditions like Alzheimer’s disease.

3 Literature Review

Novel Object Recognition (NOR) is a widely-used task to assess cognitive functions like memory. This task places mice into an arena with both a novel object and a familiar object, and their exploration is monitored. It is believed that mice with strong memory will explore the novel object more than the familiar one, a concept central to the NOR task [1]. This methodology is appealing because it does not need any external incentives or penalties and minimal habituation is required. It can be completed in a relatively short time [1]. The task procedure is divided into three phases: habituation, familiarization, and the test phase [1]. In the habituation phase, each animal freely explores the open arena with no objects present. After this process, the animal is removed from the area. In the familiarization phase, a single animal is placed in the open arena with two of the same sample objects. For a few minutes, the rodent’s movement and exploration of the objects are recorded. Finally, in the testing phase, a single animal is placed into the open area with one of the familiar objects from the familiarization phase, and one novel object. Key brain regions involved in NOR include the hippocampus and the perirhinal cortex. While the perirhinal cortex is involved in representing basic information about the familiarity or novelty of an object, the hippocampus is concerned with object memorization by encoding information about the experience of the object [1]. The parahippocampal regions of the temporal lobe, including the perirhinal, entorhinal, and inferior temporal cortices, are highly significant [3]. Notably, the perirhinal cortex’s role in object recognition decays quickly and is not sufficient for maintaining information about the object during longer retention intervals, whereas the hippocampus maintains strong novel object preference after long but not short delays [3]. NOR was first studied by Ennaceur and Delacour in 1988. They defined exploration of an object as directing the nose at a distance of less than 2 cm to the object and/or touching it with the nose, with other actions like turning around or sitting on the object not considered exploration [2]. NOR has been identified as an excellent primary screening assay for Alzheimer’s disease (AD) models, especially for evaluating spatial learning and working memory in the brain [9, 8]. Its minimal stress to mice makes it an appealing method for assessing AD’s impact on cognition. Previously, NOR has been used as a vehicle to test drugs for the treatment of Amyloid Beta and tau-impaired mice [9].

Spontaneous Location Recognition (SLR) is a behavioral assessment method used to measure spatial pattern separation ability in rodents, such as mice and rats. Unlike other pattern separation tasks, SLR varies the load on pattern separation during encoding and can be performed in standard rodent behavioral facilities, requiring only basic expertise in rodent handling [7]. The SLR task assesses the extent to which animals can discriminate and remember object locations presented during the encoding phase. By modifying the configurations of the objects in the arena and altering the distance between two of the objects, pattern separation processes can be challenged. To evaluate pattern separation performance, the amount of time the animal spends exploring the novel location relative to the familiar location is measured [7]. Impaired pattern separation has been shown to contribute to memory deficits seen in neuropsychiatric and neurodegenerative diseases, including Alzheimer’s disease. The SLR task can be used to study spatial pattern separation and memory in Alzheimer’s disease mouse models, providing valuable insights into the neurobiology of memory and cognition in various conditions, including neurodegenerative diseases [7]. Aspects such as plasticity-related mechanisms, including adult hippocampal neurogenesis, are involved in this process [1].

Any-Maze is a video-tracking and data analysis tool used in behavioral experiments, including Novel Object Recognition (NOR) tests. It allows researchers to define investigation zones around objects and determine when an animal is exploring an object. Any-Maze can provide various results, such as total time exploring the novel object, total time exploring the known object, and discrimination index. DeepLabCut is an open-source software package for markerless pose estimation of animals in videos. It uses deep learning techniques to track and analyze animal behavior, which can be applied to study exploration in NOR tests. With pose estimation of NOR and SLR video recordings, one can quantify the mice exploration of the objects. Convolutional Neural Networks (CNNs) are a class of deep learning algorithms that have shown great potential in various applications, including image classification, object detection, and semantic segmentation. They were used to analyze NOR and SLR video recordings by building off the EXPLORE open-source CNN tools [5].

In the field of behavioral experimentation, particularly in studies involving Novel Object Recognition (NOR) and Spontaneous Location Recognition (SLR) tests, technological advancements have paved the way for more precise and efficient analysis. Tools like Any-Maze offer a specialized approach for tracking and analyzing animal interactions with objects, providing insights into various behavioral patterns. Simultaneously, the open-source platform DeepLabCut (DLC) leverages the power of deep learning to enable markerless pose estimation, offering a flexible solution for researchers to quantify animal exploration. The emergence of Convolutional Neural Networks (CNNs), especially when built upon the EXPLORE open-source tools, has further expanded the horizons in video analysis, object detection, and segmentation. This paper demonstrates the efficacy of these tool in both Novel Object Recognition and Spontaneous Location Recognition whilst studying AD mice.

4 Methods

4.1 Mice

Four distinct groups are involved in this experimental design:

1. High PHF Group: Mice with increased levels of PHFs, modeling a condition analogous to advanced stages of AD.
2. Low PHF Group: Mice with lower levels of PHFs, simulating an earlier or milder stage of tau pathology.
3. B6 Control Group: A strain-specific control group using the C57BL/6 (B6) mice, providing a genetic baseline for comparison.
4. CTL (Standard Control) Group: Wild-type or unmodified mice, serving as the overall control for normal behavior.

All mice data was provided in the form of video footage in the Hussaini Lab

5 Findings

The overarching aim of this study was to develop and validate a computational pipeline for analyzing exploration behaviors in Novel Object Recognition (NOR) and Spontaneous Location Recognition (SLR) tasks. The pipeline incorporated three analytical tools: Any-Maze, DeepLabCut, and a Convolutional Neural Network (CNN) via the EXPLORE platform. Each tool served a specific purpose, and their combined application led to a rich, multidimensional analysis of the subjects' behavior.

DeepLabCut Analysis

DeepLabCut was employed for a comprehensive analysis of specific body parts of the subjects, offering an in-depth examination of their physical movements and orientations within the context of SLR and NOR video datasets. A total of 25 randomly selected frames from each video were manually labeled to identify distinct features of the mice, including the snout, tailbase, and tailtip. This labeled dataset was subsequently partitioned into training and test sets to facilitate model training.

Utilizing a ResNet (Residual Network) architecture, the training set was engaged in fine-tuning the classification of labels corresponding to specific body parts. ResNet was chosen specifically for its unique architecture that includes skip connections or shortcuts to jump over some layers. These connections enable the training of much deeper networks by alleviating the vanishing gradient

problem, leading to the enhanced performance of the network [4]. This architectural feature distinguishes ResNet from many other architectures and enables efficient training, even with a large number of layers.

The utilization of ResNet provided an advantage in capturing intricate details and patterns in the video frames, which may not have been achievable with other architectures. Its ability to train deep networks without overfitting made it well-suited for this specific application, where the fine-tuning of labels was crucial[4]. Using a successor architecture might have been superfluous for this task, as ResNet already provided the necessary depth and complexity for accurate labeling without introducing unnecessary computational burden.

Following the training, the fine-tuned model was utilized to predict labels across the remaining NOR and SLR videos, ensuring methodological uniformity and consistency throughout the analysis. Attempts were made to expand the labeling to include additional markers, such as the shoulders and hips. However, these were ultimately discarded due to the video quality not being conducive to consistent labeling. The pose estimation data generated by the above process was leveraged to construct an object recognition pipeline, enabling the creation of masks around objects of interest. The initial mask creation was facilitated using LabelMe, which allowed users to demarcate specific areas by drawing polygons. These polygons were subsequently converted into image masks.

The program was designed to quantify exploration in two distinct ways: by measuring the time spent within a range of 15 pixels around the object and by evaluating the time spent in proximity to and orientation towards the object. This approach facilitated the conversion of pose estimation data into millisecond exploration data, a more relevant and significant metric for our study.

Here are some relevant examples of DLC Analysis:

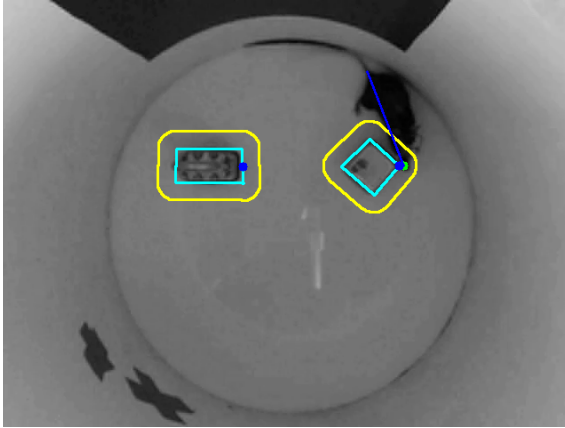


Figure 1: NOR Exploring Novel

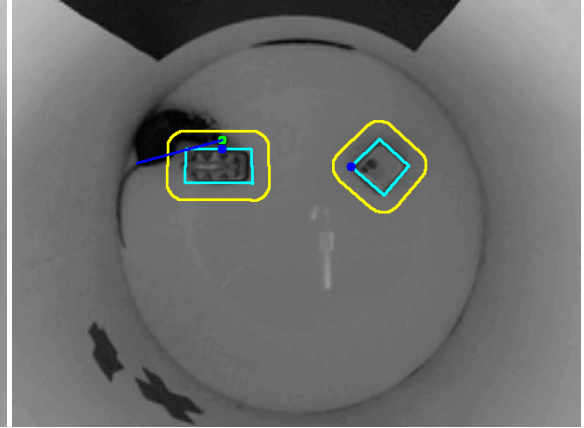


Figure 2: NOR Exploring Familiar

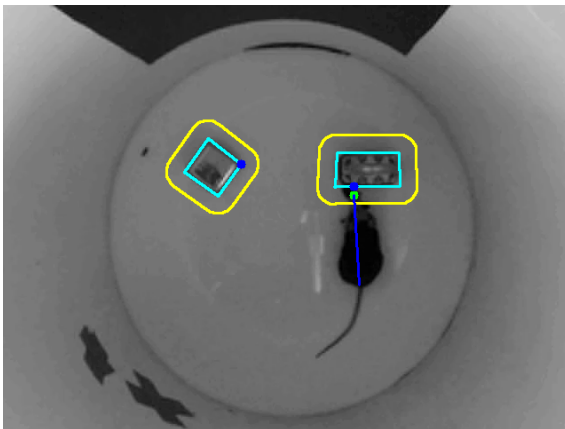


Figure 3: NOR Exploring Familiar

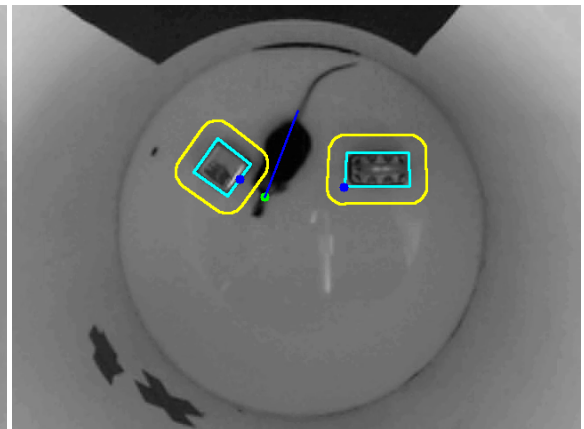


Figure 4: NOR Not Exploring

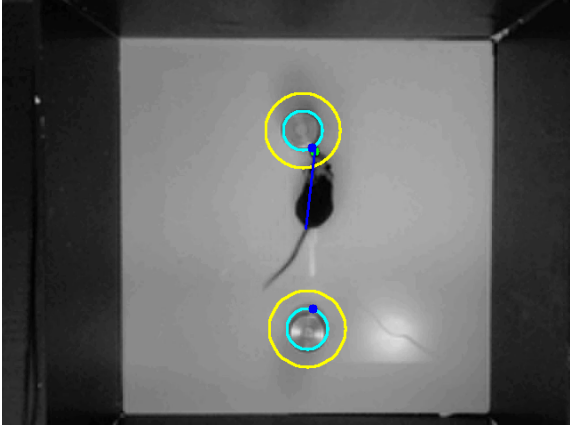


Figure 5: SLR Exploring Novel

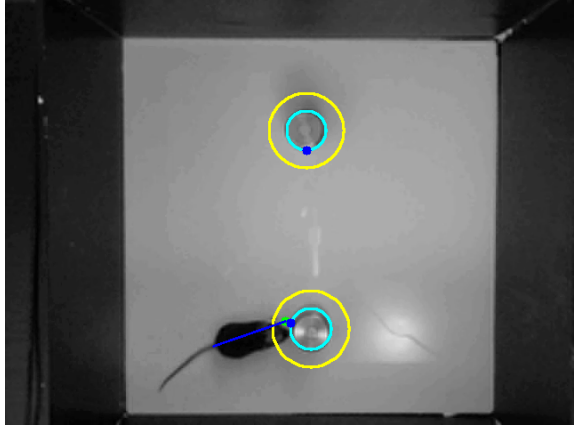


Figure 6: SLR Exploring Familiar

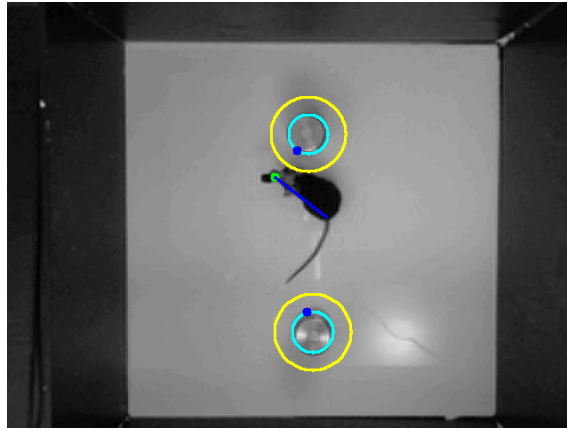


Figure 7: SLR Not Exploring

Convolutional Neural Network Analysis

EXPLORE is a deep-learning-based tool that utilizes a Convolutional Neural Network (CNN) to analyze exploration behavior, offering a nuanced understanding of the subjects' interactions with novel and familiar objects [5]. The novelty of EXPLORE lies in its flexibility, allowing users to define by themselves what object interactions of the animal are accounted for as "explorative behavior." By labeling appropriate training data, the CNN extracts features from images and classifies them into predefined classes. From a compilation of all the testing phase video data, random frames were labeled as "exploring novel," "exploring familiar," or "no," and were extracted to form a training set. Then, frames were deleted from each label to guarantee equal amounts of input data. The training and test set data were split 80:20.

For each of the three paradigms:

- NOR with a novel object on the left
- NOR with a novel object on the right
- SLR with a novel object on the top

EXPLORE was employed to train the CNN to classify frames based on the specified training data. The status of "no" refers to a state where the animal is not actively engaging in the exploration of either the novel or familiar objects, such as behaviors like resting, grooming, or other non-explorative activities. The trained CNNs were deployed to classify frames from particular videos, enabling the computation of exploration times. This comprehensive approach provided a sophisticated and flexible analysis of the subjects' detailed exploration behavior.

The models' performance can be shown below:

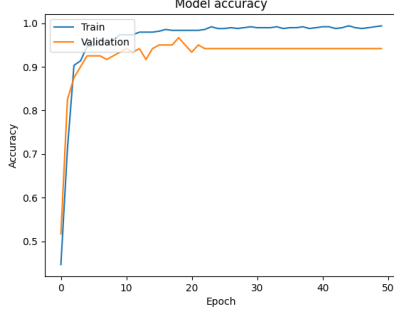


Figure 8: Right NOR Accuracy

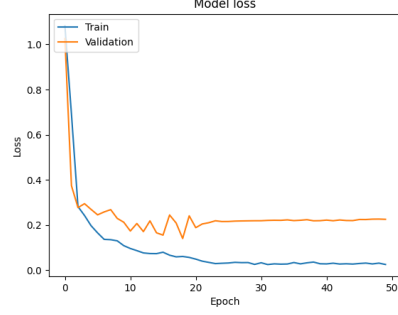


Figure 9: Right NOR Loss

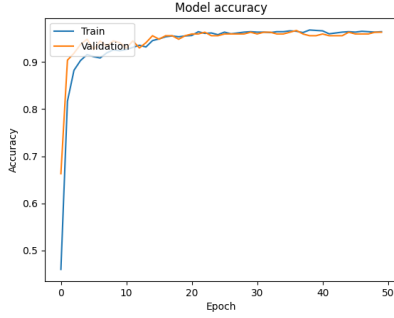


Figure 10: Left Novel Accuracy

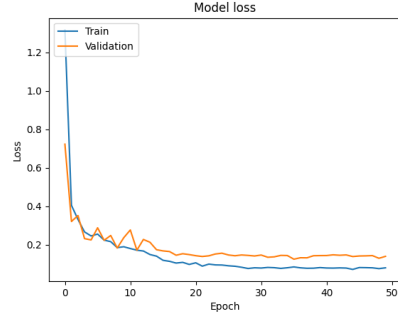


Figure 11: Left NOR Loss

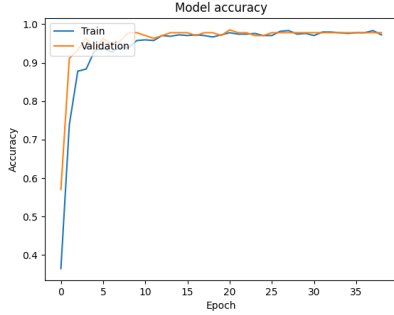


Figure 12: SLR Accuracy

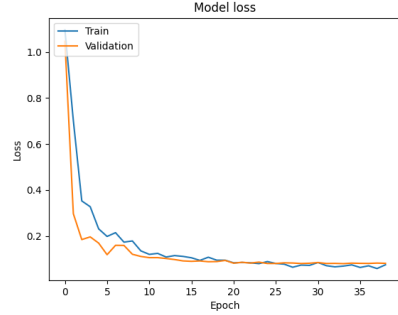


Figure 13: SLR Loss

Each model demonstrated impressive performance, achieving accuracy levels exceeding 90%. This high level of accuracy indicates the robustness and reliability of the models in classifying the given tasks. A notable exception was observed in the Right Novel training paradigm, where the testing data failed to achieve the same high levels of accuracy as the training data. This discrepancy may suggest overfitting in this specific instance or potential inconsistencies within the testing dataset. Further investigation may be required to pinpoint the underlying cause and ensure that the model generalizes well to unseen data.

5.1 Test Analysis

The study implemented various analytical techniques to investigate the exploration behaviors associated with Novel Object Recognition (NOR) and Spontaneous Location Recognition (SLR). These methodologies were designed to quantify interactions with novel and familiar objects, encompassing metrics such as time spent exploring the object and the frequency of entries and exits around it. Two principal metrics were employed to interpret the data:

1. **Discrimination Index:** This index was derived for each dataset, utilizing the equation:

$$\frac{\text{novel exploration} - \text{familiar exploration}}{\text{novel exploration} + \text{familiar exploration}}$$

The discrimination index serves as an indicator of preference for novelty over familiarity. Positive values denote a preference for novel objects, and negative values indicate a preference for familiar ones. This index provides a standardized quantitative measure to assess differential interactions with novel and familiar objects, enabling comparisons across various conditions and experiments[1].

2. **Normalized Exploration:** Alongside the discrimination index, the study also computed normalized values for both novel and familiar explorations. These normalized measurements were defined as:

$$\text{Novel Exploration Normalized} = \frac{\text{novel exploration}}{\text{novel exploration} + \text{familiar exploration}}$$

$$\text{Familiar Exploration Normalized} = \frac{\text{familiar exploration}}{\text{novel exploration} + \text{familiar exploration}}$$

These normalized values afford a relative gauge of exploration, irrespective of the total time spent. They facilitate a more detailed understanding of how subjects allocate their attention between novel and familiar objects, offering insights into cognitive mechanisms and potential biases.

Four distinct datasets were generated from the study, corresponding to various analysis techniques: anymaze, CNN, DLC (Deep Lab Cut) objects, and DLC centerpoints. Each of these datasets represents a different approach to quantifying the exploration of novel and familiar objects.

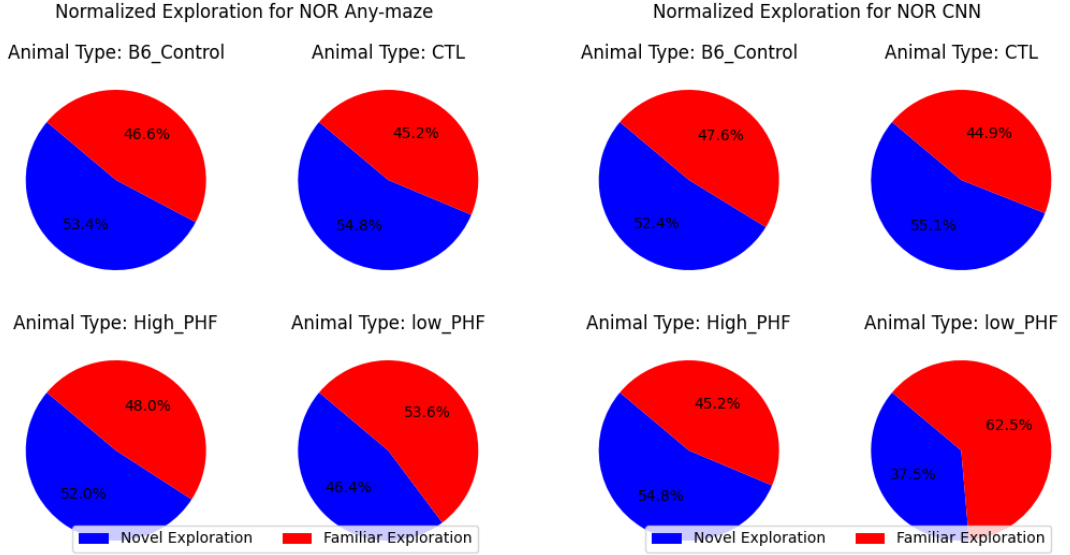


Figure 14: NOR Anymaze Pie

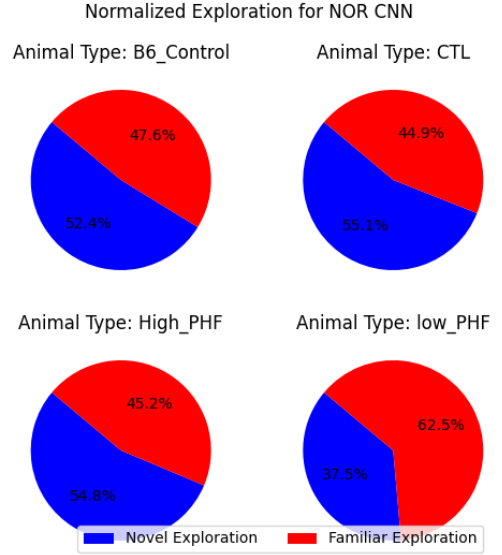


Figure 15: NOR CNN Pie

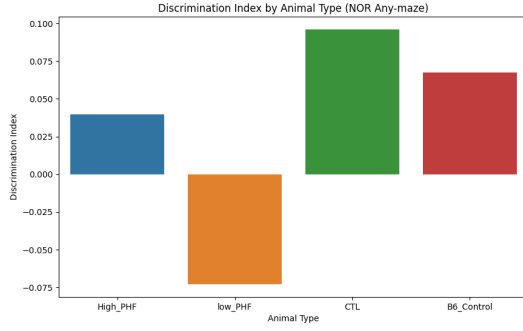


Figure 16: NOR Anymaze Plot

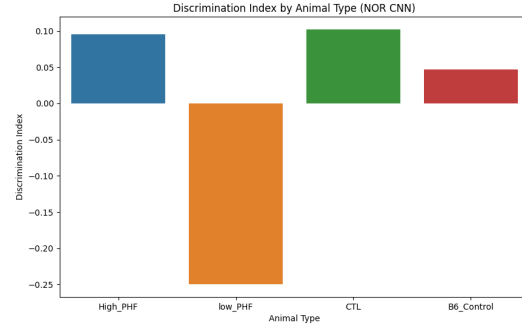


Figure 17: NOR CNN Plot

In the analysis of the two NOR samples shown, consistent behavior patterns emerged, particularly among the control groups (CTL and B6). These groups exhibited a marked preference for exploring the novel object over the familiar object. This observation is quantified by the discrimination index, a metric representing the degree to which an animal explores a novel object versus a familiar one. The index ranges between 1 and -1, with positive values indicating a preference for novel objects and negative values indicating a preference for familiar objects. However, the study also revealed an intriguing pattern within the High PHF group, which represents mice with a higher predisposition to Alzheimer's disease. Contrary to expectations, these mice also demonstrated a significant inclination to explore novel objects, closely aligning with the behavior of the control groups. This finding contrasts sharply with the behavior of the Low PHF mice, which exhibited a preference for exploring familiar objects.

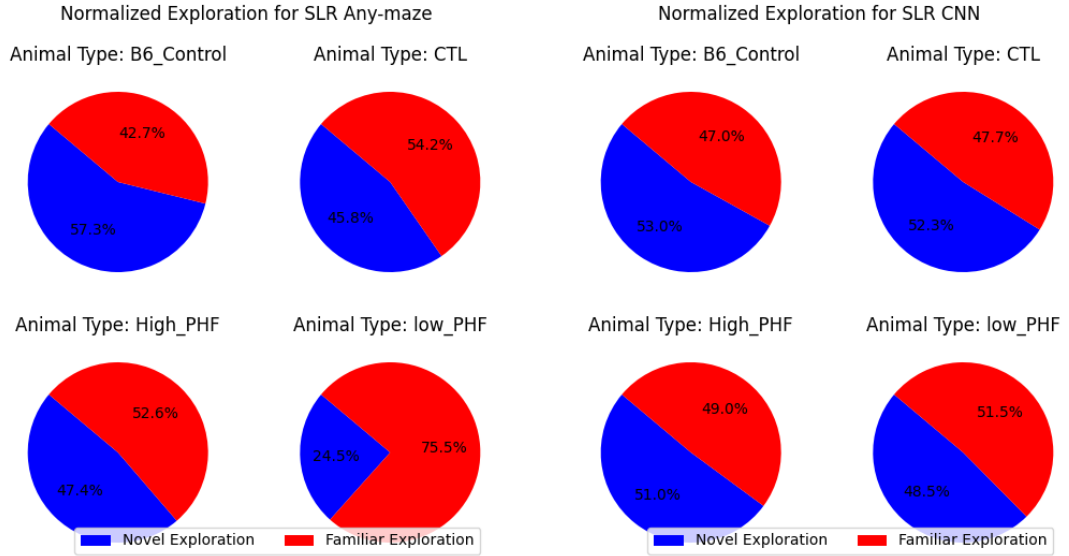


Figure 18: SLR Anymaze Pie

Figure 19: SLR CNN Pie

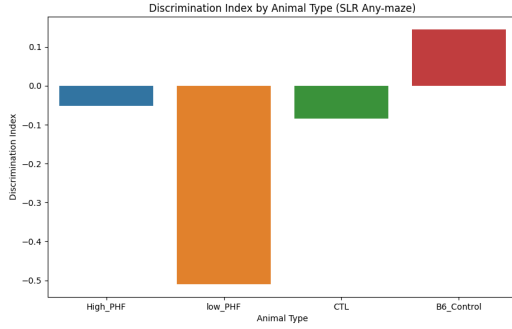


Figure 20: SLR Anymaze Plot

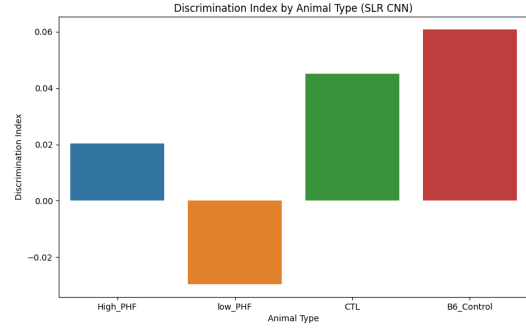


Figure 21: SLR CNN Plot

The Spontaneous Location Recognition (SLR) samples presented in this study exhibit somewhat divergent behavior patterns, which may be attributed to certain experimental errors or inconsistencies. The CNN analysis of SLR samples largely aligns with the findings from the Novel Object Recognition analysis and corroborates existing literature. Specifically, the data indicate that novel exploration is primarily observed within the control groups. However, an exception to this pattern is observed in the High PHF group, which, intriguingly, also shows a preference for novel exploration. In contrast, the anymaze data present a departure from these patterns. The control (CTL) group’s data lean slightly towards the familiar object, whereas the Low PHF group heavily leans towards familiarity. Interestingly, only the B6 control group consistently demonstrates a preference for novel exploration in the anymaze data. These discrepancies between the CNN and anymaze analyses highlight potential inconsistencies in the experimental approach or the analytical methods employed. They underscore the complexity of studying NOR and SLR behaviors, particularly in the context of Alzheimer’s disease, where nuanced and multifaceted factors may influence the outcomes.

6 Discussion

This study set out with the primary aim of establishing a robust computational pipeline for analyzing exploration behaviors in Novel Object Recognition (NOR) and Spontaneous Location Recognition (SLR). The computational pipeline successfully integrated a variety of analytical tools, including Any-Maze, DeepLabCut, and a Convolutional Neural Network through the EXPLORE platform. Each of these tools added a layer of complexity and granularity to the study, allowing for a nuanced understanding of exploration behaviors.

Among the computational challenges encountered, the limitations of DeepLabCut (DLC) were particularly significant. While DLC’s analytical framework is theoretically sound, its application was hindered by the low quality of the video recordings used for training. This issue underscores the importance of high-resolution and high-fps recordings in future iterations of this pipeline. Additionally, the non-standardization of object selection could introduce biases into the results, highlighting the need for a more computational approach to object selection in future studies.

On the statistical front, the overall analysis did not yield definitive results. However, this may be attributed to challenges like the aforementioned limitations in video quality and the potential for biases, as well as other factors like sample size and experimental conditions. Despite these setbacks, the study succeeded in implementing complex metrics like the discrimination index and normalized exploration values. These metrics were computed through the computational pipeline and offer a nuanced understanding of exploration behaviors, which will be instrumental in future research.

7 Conclusion

In conclusion, this study represents a significant stride forward in the realm of behavioral neuroscience, with a special focus on Alzheimer’s disease research. By skillfully integrating state-of-the-art computational tools such as Any-Maze, DeepLabCut, and Convolutional Neural Networks

via the EXPLORE platform, we have successfully established a multifaceted and scalable computational pipeline. This innovative pipeline not only enhances the precision and depth of behavioral analysis in Novel Object Recognition (NOR) and Spontaneous Location Recognition (SLR) tasks but also sets the stage for future, more complex studies.

A pivotal accomplishment of this research is the development of a versatile and adaptable framework. As new computational methods and tools continue to emerge, this framework can be easily modified and expanded, ensuring its ongoing relevance and efficacy in capturing the subtle complexities of animal behavior, particularly in the context of cognitive impairments like Alzheimer’s disease. This adaptability is crucial in the fast-evolving landscape of behavioral neuroscience.

Moreover, this study underlines the vital importance of high-quality data acquisition, especially in terms of video recording, for the effective application of analytical tools like DeepLabCut. It emphasizes the necessity for ongoing improvements in data collection methods and the standardization of experimental protocols, which are imperative for achieving reliable and insightful results.

The computational pipeline developed in this study lays a robust foundation for future research in this field. Its modular nature facilitates the seamless integration of additional computational tools and metrics, offering a scalable and adaptable solution to meet the evolving demands of behavioral neuroscience research. Despite facing certain challenges that impacted the statistical robustness of its findings, the study’s primary success lies in its development and validation of a comprehensive computational pipeline. This innovative approach for the behavioral analysis of NOR and SLR tasks is poised to become a benchmark for future research, offering a more nuanced and reliable methodology for data analysis in the study of neurodegenerative diseases and beyond.

8 Disclosure Of Interest

There is no competing financial interests or personal relationships that would have influenced the work reported in this paper.

9 Code Availability Statement

All code used in this research can be found at
<https://github.com/bafanaS/DLC-Object-Recognition-Analysis.git>

10 ORCID

Soham Bafana: <https://orcid.org/0009-0005-3681-1088>

References

- [1] M. Antunes and G. Biala. The novel object recognition memory: neurobiology, test procedure, and its modifications. *Cognitive Processing*, 13(2):93–110, May 2012.
- [2] A. Ennaceur and J. Delacour. A new one-trial test for neurobiological studies of memory in rats. 1: Behavioral data. *Behavioural Brain Research*, 31(1):47–59, November 1988.
- [3] R Hammond. On the delay-dependent involvement of the hippocampus in object recognition memory. *Neurobiology of Learning and Memory*, 82(1):26–34, July 2004.
- [4] Kaiming He, Xiangyu Zhang, Shaoqing Ren, and Jian Sun. Deep Residual Learning for Image Recognition. 2015.
- [5] Victor Ibañez, Laurens Bohlen, Francesca Manuella, Isabelle Mansuy, Fritjof Helmchen, and Anna-Sophia Wahl. EXPLORE: a novel deep learning-based analysis method for exploration behaviour in object recognition tests. *Scientific Reports*, 13(1):4249, March 2023.
- [6] Rodrigo Medeiros, David Baglietto-Vargas, and Frank M. LaFerla. The Role of Tau in Alzheimer’s Disease and Related Disorders: Role of Tau in AD and Related Disorders. *CNS Neuroscience & Therapeutics*, 17(5):514–524, October 2011.

- [7] Amy C. Reichelt, Cecilia P. Kramar, Olivia R. Ghosh-Swaby, Paul A. S. Sheppard, Brianne A. Kent, Pedro Bekinschtein, Lisa M. Saksida, and Timothy J. Bussey. The spontaneous location recognition task for assessing spatial pattern separation and memory across a delay in rats and mice. *Nature Protocols*, 16(12):5616–5633, December 2021.
- [8] Giulio Taglialatela, Dale Hogan, Wen-Ru Zhang, and Kelly T. Dineley. Intermediate- and long-term recognition memory deficits in Tg2576 mice are reversed with acute calcineurin inhibition. *Behavioural Brain Research*, 200(1):95–99, June 2009.
- [9] Ru Zhang, Guizhen Xue, Shaodeng Wang, Lihong Zhang, Changjie Shi, and Xin Xie. Novel Object Recognition as a Facile Behavior Test for Evaluating Drug Effects in A β PP/PS1 Alzheimer’s Disease Mouse Model. *Journal of Alzheimer’s Disease*, 31(4):801–812, September 2012.