

FineRadScore: A Radiology Report Line-by-Line Evaluation Technique Generating Corrections with Severity Scores

Alyssa Huang

Harvard University

ALYSSAHUANG@COLLEGE.HARVARD.EDU

Oishi Banerjee

Harvard University

OISHI_BANERJEE@G.HARVARD.EDU

Kay Wu

Harvard University

KAY.WU@MEDPORTAL.CA

Eduardo Pontes Reis

Stanford University

EDUARDO.REIS@EINSTEIN.BR

Hospital Israelita Albert Einstein

Pranav Rajpurkar

Harvard University

PRANAV_RAJPURKAR@HMS.HARVARD.EDU

Abstract

The current gold standard for evaluating generated chest x-ray (CXR) reports is through radiologist annotations. However, this process can be extremely time-consuming and costly, especially when evaluating large numbers of reports. In this work, we present FineRadScore, a Large Language Model (LLM)-based automated evaluation metric for generated CXR reports. Given a candidate report and a ground-truth report, FineRadScore gives the minimum number of line-by-line corrections required to go from the candidate to the ground-truth report. Additionally, FineRadScore provides an error severity rating with each correction and generates comments explaining why the correction was needed. We demonstrate that FineRadScore’s corrections and error severity scores align with radiologist opinions. We also show that, when used to judge the quality of the report as a whole, FineRadScore aligns with radiologists as well as current state-of-the-art automated CXR evaluation metrics. Finally, we analyze FineRadScore’s shortcomings to provide suggestions for future improvements. Code to run FineRadScore can be found [here](#).

1. Introduction

Artificial Intelligence (AI) is rapidly advancing in the field of medical image interpretation. In particular, models have succeeded on a range of chest X-ray (CXR) classification tasks, improving radiologists’ efficiency (Nsengiyumva et al. (2021)). Models performing the more comprehensive task of radiology report generation could further assist radiologists by describing all abnormalities in a chest X-ray in fluent text (Endo et al. (2021), Jeong et al. (2023), Nicolson et al. (2022)). However, AI models still struggle to produce reliable and accurate radiology reports. As a result, there is a need for evaluation metrics that can track progress in this space and provide detailed feedback on AI-generated reports.

Automating the evaluation of AI-generated reports is a challenging task. The current gold standard is manual evaluation by radiologists, yet evaluating a large number of re-

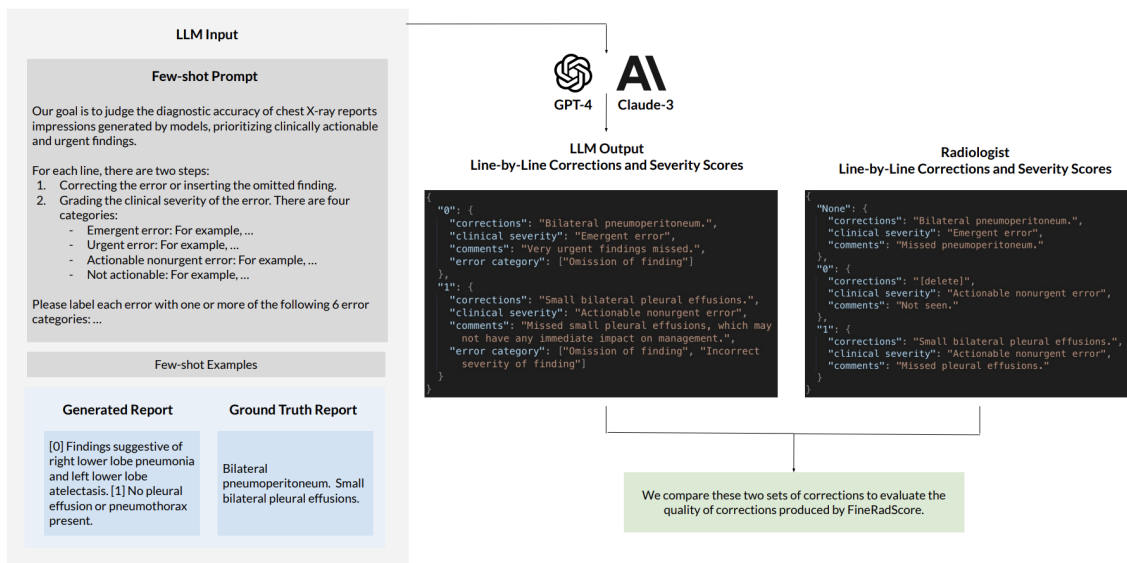


Figure 1: Overall Evaluation Framework for FineRadScore. We first prompt an LLM with a candidate report and a ground-truth report and ask for line-by-line corrections, corresponding clinical severity ratings, and comments describing each error that is corrected. Then, we feed in the pair of reports: one candidate and one ground truth. We then extract the FineRadScore-generated corrections, severity scores, and comments, and compare them with radiologist annotations. In this paper, we examine two different LLMs as backbones for the FineRadScore framework: GPT-4 and Claude-3 Opus. We will refer to FineRadScore-GPT-4 and FineRadScore-Opus to distinguish between both versions.

ports in this manner is costly and time-consuming, creating a need for automated metrics. Researchers have explored natural language generation (NLG) metrics such as BLEU (Papineni et al. (2002)) and BERTScore (Zhang et al. (2019)) to compare the generated report with a ground-truth report. However, these metrics were designed for non-medical text and fail to capture important features of medical text due to the lack of domain-specific knowledge needed to evaluate radiology reports (Yu et al. (2023), Liu et al. (2023)). For example, if a model generated a sentence such as “there is a focal lesion” and the ground-truth report stated “there is no focal lesion,” this pair would score highly on BLEU-1, even though radiologists would consider it a poor match.

As a result, recent works have developed metrics specifically designed for evaluating radiology reports, such as CheXbert vector similarity (Smit et al. (2020)), RadGraph-F1 (Jain et al. (2021)), and RadCliQ (Yu et al. (2023)). Although these metrics have shown promise, they focus on report-level evaluation, outputting a single score for the entire report. This work expands on these existing automated evaluation metrics by using large language models to offer more detailed assessments, with sentence-level feedback and fine-grained estimates of error severity.

Generalizable Insights about Machine Learning in the Context of Healthcare

In this work, we propose an evaluation framework that uses a Large Language Model (LLM) to offer sentence-level corrections, error severity ratings, and comments. We compare and contrast performance between using GPT-4 and Claude-3 Opus as our framework’s backbone. To our best knowledge, this is the first work to:

1. Develop a CXR evaluation framework that can also generate corrections
2. Develop a CXR evaluation framework that assigns clinical severity scores and error types found for each line
3. Apply and analyze Claude-3 Opus on any radiology report evaluation task.

To help evaluate our proposed framework, we also present a dataset of 100 studies with detailed sentence-level radiologist annotations. Using this new data and other existing datasets, we show that our proposed evaluation technique aligns with radiologists’ judgments, as well as identify where there is still room for improvement.

Motivation for Line-by-Line Evaluation. Combined with information on clinical error severity, the number of incorrect lines in a report helps provide a well-rounded picture of how an AI model would impact clinical practice. In particular, clinical error severity correlates with consequences on patient health, while the number of incorrect lines correlates with time needed by radiologists to correct the report. Line-by-line corrections can provide insight into which parts of a report are most flawed (e.g. “the impressions section is consistently poor-quality” or “lines mentioning fractures are frequently incorrect”). Additionally, line-by-line corrections could be used to improve models through Reinforcement Learning with AI Feedback techniques.

2. Related Work

2.1. Automated Evaluation Techniques for AI-Generated Radiology Reports.

NLG metrics such as BLEU (Papineni et al. (2002)) and BERTScore (Zhang et al. (2019)) are popular choices when evaluating AI-generated chest X-ray reports (Yu et al. (2023), Jeong et al. (2023)). BLEU measures token overlap between reports, while BERTScore compares embeddings generated by a BERT-based model that is not specific to the medical domain. In addition, other metrics specific to the medical field such as CheXbert vector similarity (Smit et al. (2020)), RadGraph-F1 (Jain et al. (2021)), and RadCliQ (Yu et al. (2023)) have been developed and been widely used to evaluate chest x-ray reports (Tu et al. (2023), Hyland et al. (2023)). CheXbert vector similarity takes in a report and produces a 14-dimensional vector, where each index indicates whether a certain pathology was found in the report, and compares cosine similarities between vectors. RadGraph-F1 extracts clinical entities and relations to create a knowledge graph from a report and computes the overlap with a ground-truth graph. RadCliQ is a weighted composite of all the aforementioned metrics (BLEU, BertScore, CheXbert, and RadGraph-F1). Among these metrics, RadCliQ has been found to be most aligned with radiologist opinions (Yu et al. (2023)).

Although these metrics can score the overall quality of a report using clinical knowledge, they are unable to give more fine-grained information. For example, they fail to explain

which parts of a report are problematic or how clinically significant an error is. Additionally, they do not indicate how the generated report can be edited to make its content align with the ground-truth report’s. This work aims to fill these gaps.

2.2. GPT-4 as an Evaluator.

Recent works have found LLMs, GPT-4 in particular, to be promising in evaluating NLG outputs. G-Eval, a GPT-4-based metric, evaluates NLG outputs using chain-of-thought prompting (Liu et al. (2023)). For a pair of candidate and ground-truth reports, G-Eval assigns a score based on the quality of the candidate report. The authors found that G-Eval is more aligned with human judgement than BLEU and BERTScore. In the medical domain, G-Rad is a similar technique, which prompts GPT-4 to identify the number of clinically significant and insignificant errors in a given report according to a set of error categories (Chaves et al. (2024)). Additionally, entailment-based prompts have been used on LLMs in the medical space to obtain a single score for each pair of reports (Zhu et al. (2023), Xie et al. (2023)).

We expand on these methods in our work by exploring whether GPT-4 can go beyond just identifying the number of errors to generating the corrections for those errors. Additionally, our work maintains the ability to detect the significance of each error by also producing a clinical severity rating for each error. Finally, we explore one of the Claude-3 models, which to the best of our knowledge, has not yet been examined in this task.

2.3. Alignment Datasets.

In order to measure how well automated metrics align with human assessments of report quality, past work has curated datasets containing radiologist scores of chest X-ray reports. Evaluation scores from automatic metrics can be compared with radiologist scores on these datasets. We introduce two such datasets, which we will use to measure how well FineRadScore aligns with human evaluations.

ReFiSco-v0. ReFiSco-v0 (Tian et al. (2023)) includes 60 chest X-ray studies, each containing several candidate reports with line-by-line expert annotations. For each line in a candidate report, radiologists were asked to either state that the line was already correct or else make a correction by rewriting or deleting the line. They were also allowed to make corrections to a report by inserting new lines. In addition, they were asked to rate each error they corrected using five clinical severity categories: ‘No error’, ‘Not actionable’, ‘Actionable nonurgent error’, ‘Urgent error’, or ‘Emergent error’ in order of increasing severity.

ReXVal. ReXVal (Yu et al. (2023)) is an expert-annotated set of 200 studies, each including an AI-generated candidate report and a radiologist-written ground-truth report. Radiologists identified errors in each candidate report by comparing it to the ground-truth report. They classified errors into six error categories: 1) False prediction of finding, 2) Omission of finding, 3) Incorrect location/position of finding, 4) Incorrect severity of finding, 5), Mention of comparison that is not present in the reference impression, 6) Omission of comparison describing a change from a previous study. When evaluating each candidate report, radiologists were asked to count the number of clinically significant and insignificant errors in each of the six error categories. In our work, when we reference these six error

categories, we will refer to them as the “ReXVal error categories” to explicitly contrast them from the clinical severity categories introduced in ReFiSco-v0.

3. Methods

3.1. FineRadScore Evaluation Framework

In Figure 1, we see an overview of FineRadScore. We feed in a pair of reports: the candidate report and a ground-truth report. Our goal is to have the LLM correctly identify the line-by-line corrections needed to make the content of the candidate report match the ground-truth report’s using the least number of changes possible. For each correction, the LLM should mark the correction as involving either the deletion, substitution, or insertion of a line. Additionally, the LLM should label the correction with one of five clinical error severity categories: “Not actionable”, “Actionable nonurgent error”, “Urgent Error”, “Emergent Error”, “Invalid Comparison”, and provide comments regarding why the correction is needed. We note that these clinical error severity categories are identical to ReFiSco-v0’s with the addition of “Invalid Comparison,” whose motivation we discuss below. To help users gain a better understanding of the types of errors being produced, we also ask it to output a list of error categories for each error it finds. We use the ReXVal error categories: “False prediction of finding”, “Omission of finding”, “Incorrect location/position of finding”, “Incorrect severity of finding”, “Mention of comparison that is not present in the reference impression”, “Omission of comparison describing a change from a previous study.” The list produced by FineRadScore may contain 0-6 error categories, as multiple error types may appear in a given line. In the body of this work, we primarily analyze how the corrections and clinical severity scores outputted by FineRadScore align with radiologists. Analyses of comments and the list of error labels according to the ReXVal categories are left to the appendix.

3.2. LLM Prompt

Our prompt describes the task, instructing the model to only include semantically relevant errors and to ignore stylistic differences. To construct a few-shot prompt, we then include five examples from our radiologist annotations, shown in Appendix A.2. We also tried using a zero-shot prompt shown in A.1, but found that it consistently underperformed the few-shot prompt. Thus, we omit zero-shot results in this paper. The full prompt is shown in Appendix A.

3.3. Intended JSON Output

We ask the LLM to output a list of corrections in JSON format shown in Appendix A.2. In practice, not every generation adhered to this format, so we would re-prompt the LLM up to a certain number of retries if the output was not formatted properly.

4. ReFiSco-v1 Dataset

The goal of the ReFiSco-v1 dataset is to ensure that FineRadScore can correct and score reports in a way that aligns with radiologists. The dataset designed is based off of the ReFiSco-v0 dataset. In this section, we describe how we created the ReFiSco-v1 dataset.

4.1. Cohort Selection

To start, we sample 100 radiology studies from the MIMIC-CXR test dataset (Johnson et al. (2019)), spanning a total of 14 subjects. We next use ClsGen (Nguyen et al. (2021)), a model with high clinical accuracy on MIMIC-CXR that is available to researchers, to produce our candidate reports.

For each study in our sample, we produce two reports: a findings report and an impressions report. First, we use ClsGen to generate the report findings. We then use GPT-4 to produce the impressions as a summarization task. We treat the findings reports and the impressions reports separately.

4.2. Annotation Process

After obtaining these 2 reports for each study, we recruited 2 radiologists to annotate subsets of the selected radiology studies. Both radiologists were given a 30-minute orientation, as well as instructions and examples to complete their annotation tasks.

The radiologists were presented with all chest X-ray images from a given study in the MIMIC-CXR dataset. Then, they were asked to annotate each report line by line. Each annotation consisted of two parts. First, the radiologist was asked to correct the line by either deleting an existing line, substituting an existing line for a new line, or inserting a new line. Then, they are asked to mark each line of the report with a clinical severity rating. In order of increasing severity, the categories are “No error”, “Invalid comparison”, “Not actionable”, “Actionable nonurgent error”, “Urgent error”, or “Emergent error”. Descriptions and examples for each of these errors were given in the instructions handout. These severity categories are the same as ReFiSco-v0’s, with the exception of “Invalid comparison,” a new category that ReFiSco-v0 grouped together with “Not actionable”. Since these models are trained on real-world datasets that make references to priors, they often make hallucinated references to priors that should be removed (Ramesh et al. (2022)). The “Invalid comparison” error category captures this failure mode. On lines where both a hallucinated prior and another error are found, we ask radiologists to rate the severity of the other error.

4.3. ReFiSco-v1 vs ReFiSco-v0

ReFiSco-v1 is similar to ReFiSco-v0 and collects radiologist annotations using the same set of instructions. The key distinctions between ReFiSco-v0 and ReFisco-v1 are that: 1) ReFiSco-v1 considers generations from the AI model ClsGen (Nguyen et al. (2021)) while ReFiSco-v0 considers generations from AI models X-REM (Jeong et al. (2023)) and CXR-RePaiR (Endo et al. (2021)) and 2) ReFiSco-v1 considers both findings and impression generations while ReFiSco-v0 only considers impressions.

5. Experiments

To assess the efficacy of FineRadScore, we treat the radiologist-corrected report version of the candidate report as the ground truth. To do this, we apply the corrections the radiologists made to the model-generated report to produce our ground-truth report. Then, we feed each pair of reports, consisting of a model-generated candidate report and its radiologist-corrected counterpart, as input into our evaluation pipeline. There were some AI-generated reports that were annotated by multiple raters. We create different pairs with each of these reports and submit each radiologist-corrected version to the LLM separately.

For each experiment, we use the evaluation framework outlined in section 3.1. Sometimes, the LLM is not able to output the correct JSON format we asked for. In those cases, we regenerate the result with up to 5 retries. In almost all of the cases, this number was sufficient to obtain a correctly-formatted output.

We decide to evaluate FineRadScore on both ReFiSco-v0 and ReFiSco-v1. ReFiSco-v0 and ReFiSco-v1 collectively contain candidate reports from four different sources (ClsGen (Nguyen et al. (2021)), X-REM (Jeong et al. (2023)), CXR-RePaiR (Endo et al. (2021)), and the MIMIC expert report (Johnson et al. (2019))). The two datasets also cover two common report sections: findings and impressions. In order to ensure that FineRadScore generalizes well across models and report types, we choose to evaluate on both of these datasets. Through our experiments, we explore whether:

1. FineRadScore correctly identifies the type of corrections (e.g. deletion, rewriting) needed for each line.
2. FineRadScore is able to correctly generate the new text that should be substituted in or inserted.
3. FineRadScore is able to determine the clinical severity rating of each error.
4. When we apply FineRadScore’s corrections to our candidate report, we end up with a report that is more aligned with the ground truth.
5. Results 1-4 hold when the candidate report is stylistically different from the ground-truth report.
6. FineRadScore performs as well as RadCliQ (Yu et al. (2023)) and G-Rad (Chaves et al. (2024)), two current state-of-the-art techniques, when scoring a report’s overall quality.

6. Results

All results shown below are generated using the few-shot prompt.

6.1. FineRadScore Correction Type Quality

We first study whether FineRadScore is able to identify the correction type needed for each line. To do so, we define categories describing the type of correction for each line in the AI-generated report: 1) No Change, 2) Delete, 3) Rewrite, or 4) Insert. The first three categories affect the original line in the generated report, while the fourth adds a new line that did not originally exist. Figure 6 in the appendix illustrates how often the radiologist performs each type of correction in ReFiSco-v1 and ReFiSco-v0.

Within each category, we measure how often FineRadScore makes the correct type of correction, matching the radiologist. For example, for each line that a radiologist chooses

not to change, we measure how often FineRadScore also makes no change. Likewise, for each line that a radiologist chooses to delete or rewrite, we measure how often FineRadScore makes the same decision on that line.

We treat insertions as a special case, as the same content can be inserted using different numbers of lines. For instance, a radiologist may insert three lines (“No focal consolidation. No pleural effusion. No pneumothorax.”), while FineRadScore may insert the exact same content using one line (“There is no focal consolidation, pleural effusion or pneumothorax”). We therefore analyze insertions at the level of whole reports, not individual lines. In other words, whenever a radiologist inserts at least one new line into a report, we measure how often FineRadScore also inserts at least one line into that report. We also measure FineRadScore’s precision while inserting lines, by examining how often FineRadScore does not insert a line when the radiologist does not insert any.

In this experiment, we feed in the following pair of reports at each time step. **Candidate Report:** Original Generated Report. **Ground-Truth Report:** Radiologist-Corrected Generated Report

Table 1: FineRadScore is able to identify the type of correction the radiologist makes at each line in the report or in the report as a whole.

Dataset	Model	Line Level Correction			Report Level Correction	
		No Change	Delete	Rewrite	Insert	No Insert
ReFiSco-v1	GPT-4	94.85%	80.60%	86.80%	68.75%	100%
	Opus	90.86%	97.39%	94.00%	93.75%	100%
ReFiSco-v0	GPT-4	91.18%	68.85%	87.07%	72.15%	100%
	Opus	73.38%	88.52%	91.81%	94.94%	100%

As shown in Table 1, FineRadScore-GPT-4 can reliably identify when a line should be left unchanged, but its performance drops when lines need to be deleted or inserted. Interestingly, in reports where a radiologist deletes a line and then makes an insertion, FineRadScore-GPT-4 often chooses to instead rewrite the line in one step, leading to low performance in identifying insertions. As a result, FineRadScore-GPT-4 may still produce a final corrected text that matches the ground truth, even though it makes different types of corrections. Overall, FineRadScore-Opus performs most consistently, achieving accuracy rates of over 90% on all correction types on ReFiSco-v1. No model inserted new lines when the radiologist did not insert any lines across both datasets, indicating robustness to this failure mode.

6.2. FineRadScore Text Rewrite and Insertion Quality

Next, we compare the lexical similarity between the new text proposed by radiologists and FineRadScore when rewriting or inserting lines. We aim to ensure that when FineRadScore decides to rewrite or insert a line, it does so in a way matching the radiologist. In Table 2, we compare the concatenation of all rewritten and inserted lines done by FineRadScore

and the concatenation of all rewritten and inserted done by the radiologist for each report. We choose to compare the concatenations because FineRadScore may choose to correct a report using different steps as the radiologist but still arrive at the same final result. For example, for a sample report ABC (original) \rightarrow ABD (corrected), FineRadScore may choose to first delete C, then insert D, while the radiologist just rewrites C as D. As shown in Table 2, FineRadScore is able to capture a majority of the text either rewritten or inserted by radiologists in its corrections.

Table 2: FineRadScore is able to capture a majority of the text either rewritten or inserted by radiologists in its corrections.

Dataset	Model	BLEU-1 Scores	BLEU-2 Scores
ReFiSco-v1	GPT-4	0.8601	0.8191
	Opus	0.8590	0.8101
ReFiSco-v0	GPT-4	0.7776	0.7398
	Opus	0.8371	0.8095

6.3. Correcting Candidate Reports Using FineRadScore

We now aim to ensure that after applying FineRadScore’s corrections to our candidate report, we obtain a corrected report that is more aligned with our ground-truth report. Therefore, we went through each correction outputted by FineRadScore and applied it to the original line in the candidate report to obtain a corrected report. We then compare the FineRadScore-corrected report and the ground-truth report in three ways. First, we compare the corrected candidate report and the ground-truth report’s text similarities using BLEU-2 scores. However, as noted in the introduction, lexical similarity need not be correlated with semantic alignment. Therefore, we use RadGraph-F1 (Jain et al. (2021)) and RadCliQ (Yu et al. (2023)), two evaluation metrics designed for chest x-ray reports and which have been shown to align well with radiologists. RadGraph-F1 measures the similarity between the sets of medically relevant words and their relationships for a pair of reports and outputs a score. The score ranges from 0 to 1, with 1 representing the highest alignment between the reports. RadCliQ is a composite metric that combines both lexical and clinical metrics for a pair of reports to output a score. The lower the score, the more aligned the pair of reports are.

To establish a baseline, we start by comparing the original uncorrected model-generated reports with the ground-truth reports using the three methods described above. Then, we compare the FineRadScore-corrected reports with the ground-truth reports to see if our corrected candidate reports are more aligned with the ground truth than the uncorrected candidate reports.

As shown in Table 3, when we apply the corrections generated by FineRadScore to our candidate report, we obtain higher BLEU scores when comparing this corrected report to our ground-truth report. This finding is consistent across models and datasets and is most prominent with FineRadScore-Opus. FineRadScore also produces corrected candidate

reports with higher RadGraph-F1 scores and lower RadCliQ scores. Therefore, according to these metrics, FineRadScore is producing higher-quality reports that are more medically aligned with the ground-truth reports.

Table 3: FineRadScore produces corrected candidate reports with higher RadGraph-F1 scores and lower RadCliQ scores. Therefore, FineRadScore is producing higher-quality reports that more aligned with the ground-truth reports.

Dataset	Model	BLEU-2 Scores		RadGraph-F1		RadCliQ	
		Base	Corrected	Base	Corrected	Base	Corrected
ReFiSco-v1	GPT-4	0.4232	0.9005	0.5376	0.9395	0.0118	-1.245
	Opus		0.9134		0.9391		-1.249
ReFiSco-v0	GPT-4	0.3994	0.8618	0.4633	0.8735	0.2285	-1.143
	Opus		0.9004		0.8935		-1.192

6.4. FineRadScore Clinical Severity Ratings

We next study whether FineRadScore is able to estimate the clinical severity of each error it corrects. Following [Tian et al. \(2023\)](#), we treat unwanted comparisons to priors as not actionable errors in quantitative analyses. When no corrections are found, we set the clinical severity rating to 0. Then, for each of the four severity ratings, we assign the numbers 1 to 4 in order of increasing severity (not actionable, actionable nonurgent error, urgent error, emergent error). We take the maximum clinical severity over all lines in a given report to produce a clinical severity score. We then compare FineRadScore’s report-level clinical severity score with the radiologist’s.

To provide a baseline for how well FineRadScore captures the severity rating for each line, we also measure how often different radiologists agree with each other, using inter-rater disagreement to provide a baseline. Inter-rater disagreement gives us a sense of the best we can expect from FineRadScore. For example, if all raters find the same error but are split between whether the error is emergent or urgent, the best we can expect from FineRadScore is to find the same error and rate it within that range of severities. In ReFiSco-v1, we did not collect annotations from different radiologists for a single report, so we cannot compute inter-rater disagreement. However, in ReFiSco-v0, for generations with multiple radiologist ratings, we choose one of the raters to be the candidate report at random and treat all other ratings as ground-truth reports. We then compare the candidate radiologist’s severity ratings to the ground truth radiologists’ to obtain inter-rater disagreement scores on ReFiSco-v0.

In [Figure 2](#), we plot the errors between radiologists’ severity ratings and FineRadScore severity ratings below for both ReFiSco-v1 and ReFiSco-v0. We also plot the inter-rater disagreement for ReFiSco-v0. Because the number of reports are different between FineRadScore and the radiologists, we divide by the number of reports.

In both graphs, there is a heavy concentration in the center for the errors found by FineRadScore, which indicates that in most cases, FineRadScore has the same severity rating

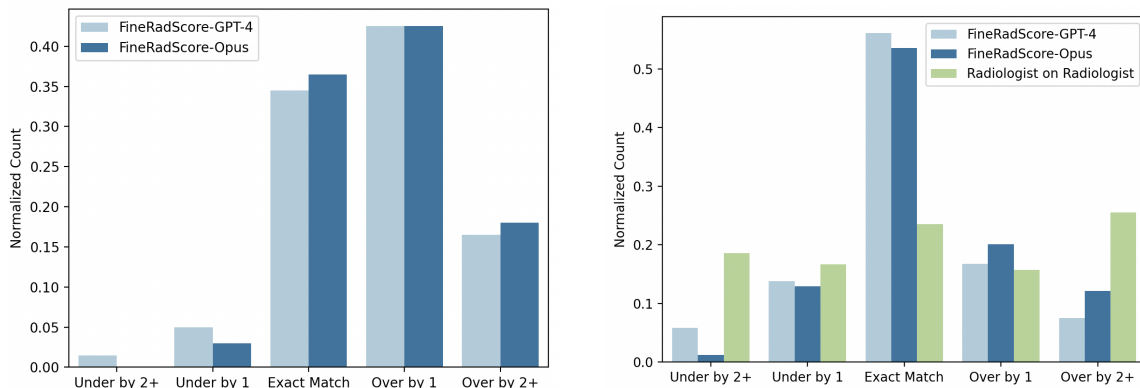


Figure 2: FineRadScore Clinical Severity Scores’ Accuracies on ReFiSco-v1 (left) and ReFiSco-v0 (right). FineRadScore’s is able to determine clinical severity scores within an error range of one category.

as the ground truth radiologist or was off by one. FineRadScore has stronger performance in rating clinical error severities in ReFiSco-v0 compared to ReFiSco-v1. We also note that FineRadScore appears to have better alignment with the ground truth radiologist compared to a candidate radiologist. However, this comparison may not be fair to the candidate radiologist since FineRadScore has access to the report corrected by the ground truth radiologist while the candidate radiologist does not. Nevertheless, these graphs indicate that FineRadScore is still able to identify the severity of each error to a significant degree. The mean absolute errors are all under 1 for FineRadScore, at around 0.84 for ReFiSco-v1 and 0.62 for ReFiSco-v0.

7. Stylistically Dissimilar Generated and Ground-truth reports

Up to this point, we have been working with candidate and ground-truth reports that are stylistically similar to each other. When moving from the candidate report to the ground-truth report, we asked radiologists to make the minimum number of corrections. Therefore, our candidate reports are stylistically similar to our ground-truth reports. However, this assumption does not always hold in practice. Therefore, we now test whether the results above hold when the candidate and ground-truth reports are stylistically dissimilar.

7.1. A Shuffled Ground Truth

We start by shuffling the order in which findings are reported in the ground truth. Then, we feed in the original candidate report with the permuted ground-truth report instead of the original ground-truth report and replicate the experiments above.

We first examine whether FineRadScore is still able to identify the type of correction needed at each line in this shuffled setting. In Table 4, we see that both FineRadScore-GPT-4 and FineRadScore-Opus are able to attain accuracies in the 80s and 90s across correction

Table 4: FineRadScore is able to identify the type of correction the radiologist makes at each line in the report or in the report as a whole even when the ground-truth report is shuffled.

Dataset	Model	Line Level Correction			Report Level Correction	
		No Change	Delete	Rewrite	Insert	No Insert
ReFiSco-v1	GPT-4	92.90%	82.46%	91.60%	80.21%	100%
	Opus	80.35%	95.15%	92.80%	95.83%	100%
ReFiSco-v0	GPT-4	86.18%	68.85%	86.64%	62.03%	100%
	Opus	61.62%	89.34%	92.24%	88.61%	100%

types on ReFiSco-v1. This trend is mostly true on ReFiSco-v0 with the exception of a few accuracies that are in the high 60s.

Table 5: FineRadScore is able to capture the majority of the text that was rewritten and inserted by the radiologist even when findings are shuffled.

Dataset	Model	BLEU-1 Scores	BLEU-2 Scores
ReFiSco-v1	GPT-4	0.8807	0.8303
	Opus	0.8223	0.7634
ReFiSco-v0	GPT-4	0.7504	0.7030
	Opus	0.8074	0.7533

Next, we examine whether the text content of the corrections produced by FineRadScore match the text of the corrections produced by the radiologist. Again, we see in Table 5 that the majority of the text that was rewritten and inserted by the radiologist is captured by FineRadScore even when findings are out of order.

Now, we apply the corrections generated by FineRadScore and see if our corrected candidate report is more aligned with the shuffled ground truth than our original candidate report. As before, we treat the comparison of the uncorrected candidate report and the ground-truth report as our baseline. In Table 6, we see a stark improvement in BLEU scores when applying the corrections generated by FineRadScore to the original candidate report. Due to the improved RadGraph-F1 and RadCliQ scores, FineRadScore’s ability to produce high quality corrections appears robust to the shuffling of sentences in the ground truth.

Finally, we examine FineRadScore’s ability to accurately rate the clinical severity of each line in the original report even with a shuffled ground-truth report. Similar as above, the mean absolute errors are all under 1, at around 0.89 for ReFiSco-v1 and 0.67 for ReFiSco-v0.

Table 6: FineRadScore is able to produce a corrected report that aligns well with the ground-truth report even when findings are shuffled.

Dataset	Model	BLEU-2 Scores		RadGraph-F1		RadCliQ	
		Base	Corrected	Base	Corrected	Base	Corrected
ReFiSco-v1	GPT-4	0.4043	0.8703	0.5337	0.9363	0.0711	-1.162
	Opus		0.8601		0.9238		-1.1332
ReFiSco-v0	GPT-4	0.3824	0.8419	0.4533	0.8610	0.2720	-1.068
	Opus		0.8627		0.8680		-1.096

7.2. A Paraphrased Candidate Report

Now, we want to see how FineRadScore performs when the candidate and ground-truth reports are stylistically distinct. Therefore, we ask GPT-4 to rephrase a subset of report generations in ReFiSco-v1 and had radiologists double-check that the paraphrases maintain the same semantic meaning. We note that in our few-shot prompt, we emphasize to FineRadScore to ignore changes in phrasing and focus on semantic differences that may lead to clinically relevant errors.

As before, we start by examining if FineRadScore is able to identify the type of correction needed at each line when we paraphrase the candidate report. In Table 7, we see that FineRadScore struggles more heavily in the “No Change” category. For both FineRadScore-Sonnet and FineRadScore-Opus, performance in the “No Change” category is quite low. On the other hand, the “Rewrite” category has the highest accuracy across models. We conclude that in the paraphrased setting, it seems likely that FineRadScore is prone to over-rewriting. In particular, lines which should not have been rewritten may be getting rewritten due to stylistic differences, rather than differences in clinical meaning.

Table 7: FineRadScore falls short in identifying lines that should be left alone. It appears to be prone to over-rewriting.

Dataset	Model	Line Level Correction			Report Level Correction	
		No Change	Delete	Rewrite	Insert	No Insert
ReFiSco-v1	GPT-4	72.34%	62.96%	87.10%	78.57%	100%
	Opus	51.06%	92.59%	93.55%	100.00%	100%

Similarly, when comparing the BLEU scores in Table 8, FineRadScore has lower rewrite and insertion BLEU scores in the paraphrased setting. However, when looking at the lines that FineRadScore rewrites, there are many instances where the rewrites are stylistically different but semantically similar to each other. For example, the radiologist rewrites some line to “Stable post CABG with intact sternal wires” and FineRadScore rewrites the same line as “Sternal wires from coronary artery bypass graft surgery are stable”. As a result,

the BLEU score for this rewritten line may be low for FineRadScore, but it is still able to capture the meaning of the correction.

Table 8: FineRadScore has lower BLEU scores when comparing the concatenation of rewritten and inserted text to the radiologist.

Dataset	Model	BLEU-1 Score	BLEU-2 Scores
ReFiSco-v1	GPT-4	0.7060	0.6640
	Opus	0.6224	0.5847

We see in Table 9 that all models result in corrected candidate reports with better RadGraph-F1 and RadCliQ scores. These results are positive indications that the corrections generated by FineRadScore still produce a high quality corrected report, even in a paraphrased setting. However, they may be prone to producing extraneous corrections based on stylistically differences that are not clinically relevant.

Table 9: All models result in corrected candidate reports with better RadGraph-F1 and RadCliQ scores in the paraphrased setting.

Dataset	Model	BLEU-2 Scores		RadGraph-F1		RadCliQ	
		Base	Corrected	Base	Corrected	Base	Corrected
ReFiSco-v1	GPT-4	0.0705	0.4753	0.2252	0.6858	0.9811	-0.4583
	Opus		0.6150		0.7587		-0.7417

Finally, we analyze FineRadScore’s ability to identify the clinical severity of the error associated with each line (or if no error exists). The mean absolute errors are all under 1 for FineRadScore, at around 0.89 for GPT-4 and 0.58 for Opus. FineRadScore is still able to identify the clinical error severity associated with each line within one category even when the candidate and ground-truth reports are stylistically different.

8. Alignment with Radiologists on ReXVal

RadCliQ (Yu et al. (2023)), G-Rad (Chaves et al. (2024)), and GREEN (Ostmeier et al. (2024)) are some current state-of-the-art methodologies for chest x-ray report evaluation and have been shown to be more aligned with radiologist evaluation than prior metrics. RadCliQ is a composite metric that outputs a single score based off of lexical and clinical metrics for each pair of reports. Given a candidate and ground truth, G-Rad prompts GPT-4 to output the number of errors in each of the six ReXVal error categories found in a candidate report. Additionally, G-Rad instructs GPT-4 to distinguish between clinically significant and clinically insignificant errors, echoing the instructions given to the radiologists while creating the ReXVal dataset. Similarly, given a candidate and ground truth report pair, GREEN outputs a score based on the number of clinically significant errors in each of the

six ReXVal error categories. All of these metrics were benchmarked on the ReXVal dataset (Yu et al. (2023)) in order to compare alignment with radiologists with prior evaluation techniques. We therefore also compare FineRadScore’s alignment with radiologists using ReXVal.

However, one key distinction between FineRadScore and prior evaluation metrics is that FineRadScore provides line-level, not report-level, scores for each candidate report. Therefore, we need a way to get a single score for each candidate report from the corrections. Just as in our experiments, when no corrections are found, we set the clinical severity rating to 0. Then, for each of the four severity ratings, we assign the numbers 1 to 4 in order of increasing severity (not actionable, actionable nonurgent error, urgent error, emergent error). We then sum the clinical severity scores of all corrections generated for a given report to get a single FineRadScore score for each report. ReXVal contains radiologist annotations counting the total number of errors in each of six error categories for each report pair. These are the radiologist scores for each report. Using the same bootstrapping technique with 1000 resamples as Yu et al. (2023), we measure the Kendall tau correlation between FineRadScore’s scores and the radiologist scores.

The ReXVal dataset contains radiologist annotations for 200 pairs of candidate and ground-truth reports. Since RadCliQ used part of the ReXVal dataset to train their composite metric, they evaluated on a held out test set of 40 annotated pairs. Therefore, we follow their methodology and compute the Kendall tau correlation on that same held out test set. The Kendall tau correlations, as well as the 95% confidence intervals are reported in Table 10. BLEU, BERTScore, CheXbert vector similarity, RadGraph, and RadCliQ numbers are taken from Yu et al. (2023). G-Rad evaluates their methodology on the full ReXVal dataset. Therefore, we use the same technique to obtain Kendall tau scores on the full ReXVal dataset and report the results in Table 11.

Table 10: FineRadScore has a Kendall tau b correlation comparable to that of RadCliQ and other baseline metrics on ReXVal’s held-out test set of 40 data points.

	Kendall tau b correlation
BLEU	0.414 (95% CI, 0.156 0.635)
BERTScore	0.505 (95% CI, 0.273 0.671)
CheXbert	0.537 (95% CI, 0.330 0.717)
RadGraph	0.528 (95% CI, 0.357 0.687)
RadCliQ	0.615 (95% CI, 0.450 0.749)
FineRadScore (GPT-4)	0.701 (95% CI, 0.523 0.841)
FineRadScore (Claude-3 Opus)	0.737 (95% CI 0.593 0.850)

Table 10 shows that FineRadScore is just as aligned (if not slightly more aligned) with radiologists as RadCliQ on this held out test set. Table 11 shows that FineRadScore is slightly less aligned as G-Rad on the ReXVal dataset. However, in both of these comparisons, the 95% CI between FineRadScore and RadCliQ overlap, indicating that FineRadScore is not

Table 11: FineRadScore has a Kendall tau b correlation lower, but still with overlapping 95% CI, compared to that of G-Rad and summed error counts on GREEN on the full ReXVal dataset. FineRadScore has a higher Kendall tau b correlation compared to the GREEN score.

	Kendall tau b correlation
G-Rad	0.76 (95% CI, 0.70 0.80)
GREEN	0.63 (95% CI, 0.56 0.69)
Error counts GREEN	0.79 (95% CI, 0.74 0.83)
FineRadScore (GPT-4)	0.700 (95% CI, 0.631 0.756)
FineRadScore (Claude-3 Opus)	0.738 (95% CI, 0.680 0.788)

significantly underperforming on ReXVal compared to existing state of the art evaluation techniques.

The key advantage that FineRadScore brings compared to these existing techniques is its ability to offer much more fine-grained information. Both RadCliQ and G-Rad offer a single score for a candidate report. GREEN also outputs the most representative error explanations for each of the six ReXVal categories. However, FineRadScore goes one step further by giving us the line-by-line corrections, as well as clinical severity ratings. As a result, FineRadScore can help us narrow down which portion of the report is most problematic, as well as what steps can be taken to correct the errors in the candidate report.

9. Qualitative Analysis of FineRadScore

We experimented with FineRadScore when the candidate and ground-truth reports are semantically identical. When the candidate and ground truths are identical, both FineRadScore-GPT-4 and FineRadScore-Opus output no corrections. When they are shuffled versions of each other, FineRadScore-GPT-4 outputs 0.03 corrections per report and FineRadScore-Opus 0.315 per report. When they are paraphrased versions of each other, FineRadScore-GPT-4 outputs 0.75 corrections per report and FineRadScore-Opus 3.5 per report. The majority of these corrections were labeled as non-actionable. However, these results indicate how FineRadScore is still susceptible to generating corrections purely due to changes in phrasing.

10. Limitations

Future work may focus on improving performance when the ground truth differs significantly stylistically from the candidate report. As shown in Section 9, FineRadScore is still prone to identifying corrections based off of syntactic differences instead of clinically significant differences. Even so, it is able to achieve comparable performance to state of the art evaluation techniques such as RadCliQ (Yu et al. (2023)) and G-Rad (Chaves et al. (2024)). Therefore, we recommend future work focus on finding a method that can solely identify clinically relevant errors, while ignoring differences in phrasing. We hypothesize that such

a finding would improve the performance of LLM-based evaluation techniques across the board.

Additionally, we recommend creating a new dataset using the expert report from the MIMIC-CXR dataset as ground truth. In other words, for a set of model generations, expert annotations regarding the minimum number of corrections needed to transform the model-generated report to be semantically equivalent to the MIMIC-CXR ground-truth report should be collected. In this way, we would be able to evaluate FineRadScore, as well as other evaluation techniques, on this dataset. Evaluation techniques performing strongly on this dataset would likely be most useful to researchers who could now evaluate their own model generations against the MIMIC-CXR ground-truth report.

Ethical Considerations In all of our experiments, we used Azure’s OpenAI endpoints to access GPT-4. We used Anthropic’s endpoints to access Claude-3 Opus. We also signed and adhered to the PhysioNet Data Use Agreement to ensure the confidentiality of the MIMIC-CXR data.

11. Conclusion

In this work, we introduce FineRadScore, a LLM-based method of obtaining fine-grained evaluation of AI-generated chest X-ray reports. We demonstrate that FineRadScore is able to generate corrections in a way that aligns with radiologists and show how it can be used to evaluate generations in practice. We also demonstrate that FineRadScore’s report-level scores align with radiologists’ approximately as well as RadCliQ and G-Rad, two state of the art evaluation techniques. However, more work is needed to evaluate pairs of reports that are stylistically dissimilar, in addition to having different content.

References

- Juan Manuel Zambrano Chaves, Shih-Cheng Huang, Yanbo Xu, Hanwen Xu, Naoto Usuyama, Sheng Zhang, Fei Wang, Yujia Xie, Mahmoud Khademi, Ziyi Yang, Hany Awadalla, Julia Gong, Houdong Hu, Jianwei Yang, Chunyuan Li, Jianfeng Gao, Yu Gu, Cliff Wong, Mu Wei, Tristan Naumann, Muhao Chen, Matthew P. Lungren, Serena Yeung-Levy, Curtis P. Langlotz, Sheng Wang, and Hoifung Poon. Training small multi-modal models to bridge biomedical competency gap: A case study in radiology imaging, 2024.
- Mark Endo, Rayan Krishnan, Viswesh Krishna, Andrew Y. Ng, and Pranav Rajpurkar. Retrieval-based chest x-ray report generation using a pre-trained contrastive language-image model. In *Machine Learning for Health (ML4H)*, pages 209–219, 2021.
- Stephanie L. Hyland, Shruthi Bannur, Kenza Bouzid, Daniel C. Castro, Mercy Ranjit, Anton Schwaighofer, Fernando Pérez-García, Valentina Salvatelli, Shaury Srivastav, Anja Thieme, Noel Codella, Matthew P. Lungren, Maria Teodora Wetscherek, Ozan Oktay, and Javier Alvarez-Valle. Maira-1: A specialized large multimodal model for radiology report generation. 2023. URL <https://arxiv.org/pdf/2311.13668.pdf>.

- Saahil Jain, Ashwin Agrawal, Adriel Saporta, Steven QH Truong, Du Nguyen Duong, Tan Bui, Pierre Chambon, Yuhao Zhang, Matthew P. Lungren, Andrew Y. Ng, Curtis P. Langlotz, and Pranav Rajpurkar. Radgraph: Extracting clinical entities and relations from radiology reports. 2021. URL [arXiv:2106.14463](https://arxiv.org/abs/2106.14463).
- Jaehwan Jeong, Katherine Tian, Andrew Li, Sina Hartung, Subathra Adithan, Fardad Behzadi, Juan Calle, David Osayande, Michael Pohlen, and Pranav Rajpurkar. Multimodal image-text matching improves retrieval-based chest x-ray report generation. In Ipek Oguz, Jack Noble, Xiaoxiao Li, Martin Styner, Christian Baumgartner, Mirabela Rusu, Tobias Heinmann, Despina Kontos, Bennett Landman, and Benoit Dawant, editors, *Medical Imaging with Deep Learning*, volume 227 of *Proceedings of Machine Learning Research*, pages 978–990. PMLR, 10–12 Jul 2023. URL <https://proceedings.mlr.press/v227/jeong24a.html>.
- Alistair E. W. Johnson, Tom J. Pollard, Nathaniel R. Greenbaum, Matthew P. Lungren, Chih ying Deng, Yifan Peng, Zhiyong Lu, Roger G. Mark, Seth J. Berkowitz, and Steven Horng. Mimic-cxr-jpg, a large publicly available database of labeled chest radiographs. 2019. URL <https://arxiv.org/abs/1901.07042>.
- Yang Liu, Dan Iter, Yichong Xu, Shuohang Wang, Ruochen Xu, and Chenguang Zhu. G-eval: Nlg evaluation using gpt-4 with better human alignment. 2023. URL <https://arxiv.org/abs/2303.16634>.
- Hoang Nguyen, Dong Nie, Taivanbat Badamdorj, Yujie Liu, Yingying Zhu, Jason Truong, and Li Cheng. Automated generation of accurate & fluent medical X-ray reports. In *Proceedings of the 2021 Conference on Empirical Methods in Natural Language Processing*, pages 3552–3569, Online and Punta Cana, Dominican Republic, November 2021. Association for Computational Linguistics. doi: 10.18653/v1/2021.emnlp-main.288. URL <https://aclanthology.org/2021.emnlp-main.288>.
- Aaron Nicolson, Jason Dowling, and Bevan Koopman. Improving chest x-ray report generation by leveraging warm starting. 2022. URL <https://doi.org/10.48550/arXiv.2201.09405>.
- Ntwali Placide Nsengiyumva, Hamidah Hussain, Olivia Oxlade, Arman Majidulla, Ahsana Nazish, Dick Menzies Aamir J Khan, Faiz Ahmad Khan, and Kevin Schwartzman. Triage of persons with tuberculosis symptoms using artificial intelligence–based chest radiograph interpretation: A cost-effectiveness analysis. In *Open Forum Infectious Diseases*, 2021. URL <https://doi.org/10.1093/ofid/ofab567>.
- Sophie Ostmeier, Justin Xu, Zhihong Chen, Maya Varma, Louis Blankemeier, Christian Bluethgen, Arne Edward Michalson, Michael Moseley, Curtis Langlotz, Akshay S Chaudhari, and Jean-Benoit Delbrouck. Green: Generative radiology report evaluation and error notation, 2024.
- Kishore Papineni, Salim Roukos, Todd Ward, and Wei-Jing Zhu. Bleu: a method for automatic evaluation of machine translation. pages 311–318, 2002. URL <https://doi.org/10.3115/1073083.1073135>.

- Vignav Ramesh, Nathan A. Chi, and Pranav Rajpurkar. Improving radiology report generation systems by removing hallucinated references to non-existent priors. In Antonio Parziale, Monica Agrawal, Shalmali Joshi, Irene Y. Chen, Shengpu Tang, Luis Oala, and Adarsh Subbaswamy, editors, *Proceedings of the 2nd Machine Learning for Health symposium*, volume 193 of *Proceedings of Machine Learning Research*, pages 456–473. PMLR, 28 Nov 2022. URL <https://proceedings.mlr.press/v193/ramesh22a.html>.
- Akshay Smit, Saahil Jain, Pranav Rajpurkar, Anuj Pareek, Andrew Ng, and Matthew Lungren. Combining automatic labelers and expert annotations for accurate radiology report labeling using BERT. In Bonnie Webber, Trevor Cohn, Yulan He, and Yang Liu, editors, *Proceedings of the 2020 Conference on Empirical Methods in Natural Language Processing (EMNLP)*, pages 1500–1519, Online, November 2020. Association for Computational Linguistics. doi: 10.18653/v1/2020.emnlp-main.117. URL <https://aclanthology.org/2020.emnlp-main.117>.
- Katherine Tian, Sina J Hartung, Andrew A Li, Jaehwan Jeong, Fardad Behzadi, Juan Calle-Toro, Subathra Adithan, Michael Pohlen, David Osayande, and Pranav Rajpurkar. Refisco: Report fix and score dataset for radiology report generation (version 0.0). 2023. URL <https://doi.org/10.13026/cneg-zk64>.
- Tao Tu, Shekoofeh Azizi, Danny Driess, Mike Schaekermann, Mohamed Amin, Pi-Chuan Chang, Andrew Carroll, Chuck Lau, Ryutaro Tanno, Ira Ktena, Basil Mustafa, Aakanksha Chowdhery, Yun Liu, Simon Kornblith, David Fleet, Philip Mansfield, Sushant Prakash, Renee Wong, Sunny Virmani, Christopher Semturs, S Sara Mahdavi, Bradley Green, Ewa Dominowska, Blaise Aguera y Arcas, Joelle Barral, Dale Webster, Greg S. Corrado, Yossi Matias, Karan Singhal, Pete Florence, Alan Karthikesalingam, and Vivek Natarajan. Towards generalist biomedical ai. 2023. URL <https://arxiv.org/pdf/2307.14334.pdf>.
- Yiqing Xie, Sheng Zhang, Zelalem Gero Hao Cheng, Cliff Wong, Tristan Naumann, and Hoifung Poon. Enhancing medical text evaluation with gpt-4. 2023. URL <https://arxiv.org/pdf/2311.09581.pdf>.
- Feiyang Yu, Mark Endo, Rayan Krishnan, Ian Pan, Andy Tsai, Eduardo Pontes Reis, Eduardo Kaiser Ururahy Nunes Fonseca, Henrique Min Ho Lee, Zahra Shakeri Hossein Abad, Andrew Y. Ng, Curtis P. Langlotz, Vasantha Kumar Venugopal, and Pranav Rajpurkar. Evaluating progress in automatic chest x-ray radiology report generation. 2023. URL <https://doi.org/10.1016/j.patter.2023.100802>.
- Tianyi Zhang, Felix Wu Varsha Kishore, Kilian Q. Weinberger, and Yoav Artzi. Bertscore: Evaluating text generation with bert. 2019. URL <https://arxiv.org/abs/1904.09675>.
- Qingqing Zhu, Xiuying Chen, Qiao Jin, Benjamin Hou, Tejas Sudharshan Mathai, Pritam Mukherjee, Xin Gao, Ronald M. Summers, and Zhiyong Lu. Leveraging professional radiologists’ expertise to enhance llms’ evaluation for radiology reports. 2023. URL <https://arxiv.org/pdf/2303.16634.pdf>.

Appendix A.

Appendix A. Prompt Instructions

We started off each prompt with the same set of prompt instructions, which is a redacted version of the instructions we gave the radiologists. Both few-shot and zero-shot prompts began with these set of instructions. These instructions are shown below.

Instructions: Return only a json object for this radiology report, with a key-value pair for every line.

Each line starts with a numerical id.

The key will be the id.

The value will be another JSON object.

Our goal is to judge prioritizing clinically actionable and urgent findings. We are looking to determine the accuracy of report impressions generated by models or reports from human radiologists.

For each line, there are two steps:

1. Correcting the error or inserting the omitted finding:
 - To make a change to an existing line, copy and paste the text and make your changes.
 - To add a new line, simply insert a row as necessary.
 - To entirely delete a line, please enter [delete].
2. Grading the clinical severity of the error.

There are four categories:

For determining the error urgency, please take into consideration that the patients were seen in an ICU setting.

- Emergent error: Findings that suggest a need for immediate intervention, with significant impact on patient's health
Examples: missed or incorrectly called tension pneumothorax, pneumoperitoneum, significantly malpositioned endotracheal tube
- Urgent error: Findings where failure to act may adversely affect patient health, that require urgent (but not immediate attention) and that if not acted on, may worsen over time and likely result in an adverse patient outcome
Examples: missed pneumonia, mildly malpositioned endotracheal tube, incorrect anatomy for pleural effusion
- Actionable nonurgent error: Findings that likely do not require action in the short term but if not acted upon may reasonably impact the patient's health.
Examples: pulmonary nodule which requires follow up, presence of emphysematous changes
- Not actionable: Findings that reasonably would not likely have an impact on the patient's management. Errors that contain a reference to a nonexistent

comparison but otherwise contain the correct information are in this category.

Examples: chronic appearing rib fracture, describing a consolidation as stable when there is no comparison, description of pulmonary edema as asymmetric rather than bilateral (as this would not change management)

Only include semantically relevant errors. For example, if the report states "no pleural effusion" and the model generates "no pleural effusion is seen", this is not an error.

However, incorrect sizes or locations of findings are considered errors. Incorrect findings should be removed and omitted findings should be added.

Please label each error with one or more of the following 6 error categories:

1. False prediction of finding
2. Omission of finding
3. Incorrect location/position of finding
4. Incorrect severity of finding
5. Mention of comparison that is not present in the reference impression
6. Omission of comparison describing a change from a previous study

Please format your output as a JSON object as shown below.

A.1. Zero-Shot Prompt

In our zero-shot prompt, we simply provide some rules regarding how we want our output to be formatted, so that we can postprocess the GPT-4 generations. The zero-shot prompt is as follows.

Results should be a JSON object in the format of:

```
{
  "1": {
    "corrections": "Corrected sentence 1 from report.",
    "clinical severity": "severity of error",
    "comments": "any comments regarding the correction",
    "error category": ["list of error categories"]
  },
  "5": {
    "corrections": "Corrected sentence 5 from report.",
    "clinical severity": "severity of error",
    "comments": "any comments regarding the correction",
    "error category": ["list of error categories"]
  },
  "None": {
    "corrections": "Insertion of new sentence for report.",
```

```

    "clinical severity": "severity of error",
    "comments": "any comments regarding the correction",
    "error category": ["list of error categories"]
  },
}

```

A.2. Few-Shot Prompt

Below we include the five shots used as part of our few-shot prompt. One of the shots contained no corrections and was designed to emphasize the fact that stylistic changes should be ignored. Three of the examples were taken from the sample examples given to all radiologists at the start of the annotation process. One of the example was taken from one of the ReFiSco-v1 annotations that was excluded from the dataset in all experiments. The examples are shown in Figure 3.

Impression lines	Corrected line	Error clinical severity	Comments
___M with trauma, evaluate for injuries, pneumothorax. Left lower lung platelike atelectasis.	Right lower lung consolidation, either pneumonia, aspiration, or possibly pulmonary contusions from recent trauma.	Actionable nonurgent error	Incorrect anatomy, which although incorrect likely does not result in urgent change in management
Cardiomegaly	[delete]	Actionable nonurgent error	Cardiomegaly not present, which may result in unnecessary work up but likely not urgent in nature
	No evidence of displaced rib fracture or pneumothorax.	Not actionable	Given the indication, this was added.
Three left lung nodules concerning for metastatic disease.	Two left lung nodules concerning for metastatic disease.	Actionable nonurgent error	Incorrect number of nodules, may affect longer term follow up
Multiple lung nodules	Left basilar opacity could represent atelectasis or consolidation.	Urgent error	In the context of neutropenic fever, missed possible pneumonia is an urgent finding
	[delete]	Not actionable	Repetitive
Stable position of endotracheal tube projects 2.2 cm above the carina	Endotracheal tube projects approximately 2.2 cm above the carina.	Not actionable	Report references a non existent comparison (semantic error) but this does not affect management
Minimal atelectasis at the right lung base.	Minimal atelectasis at the left and right lung bases.	Not actionable	Incorrect anatomy of atelectasis would not affect management
Moderate cardiomegaly.			
Pulmonary edema	[delete]	Urgent error	Incorrect call of pulmonary edema is a potentially urgent error if diuretics are started, for example
The presence of a minimal left pleural effusion cannot be excluded.			
Endotracheal tube is in standard position. A nasogastric tube is seen coursing into the stomach with tip in the stomach.	[delete]	Actionable nonurgent error	no gastric tube
Heart size is normal.			
Lungs are clear.			
No pleural effusion or pneumothorax.			
Enteric tube tip is in the stomach.	[delete]	Actionable nonurgent error	no enteric tube

Figure 3: Few-shot Examples Used by FineRadScore

The few-shot examples look as follows:

Example 1:

```

Generated Text: [0] ___M with trauma, evaluate for injuries, pneumothorax.
[1] Left lower
lung platelike atelectasis.
[2] Cardiomegaly.

```

Ground Truth Text: Right lower lung consolidation, either pneumonia,

FINERADSCORE

aspiration, or possibly pulmonary contusions from recent trauma.
No evidence of displaced rib fracture or pneumothorax.

Corrections:

```
{
  "0": {
    "corrections": "Right lower lung consolidation, either pneumonia,
                  aspiration, or
                  possibly pulmonary contusions from recent trauma.",
    "clinical severity": "Actionable nonurgent error",
    "comments": "Incorrect anatomy, which although incorrect
                likely does
                not result in urgent change in management",
    "error category": ["Incorrect location/position of finding"]
  },
  "2": {
    "corrections": "[delete]",
    "clinical severity": "Actionable nonurgent error",
    "comments": "Cardiomegaly not present, which may result in
                unnecessary
                work up but likely not urgent in nature",
    "error category": ["False prediction of finding"]
  },
  "None": {
    "corrections": "No evidence of displaced rib fracture
                  or pneumothorax.",
    "clinical severity": "Not actionable",
    "comments": "Given the indication, this was added.",
    "error category": ["Omission of finding"]
  }
}
```

Example 2:

Generated Text:

```
[0] Three left lung nodules concerning for metastatic disease.
[1] Multiple lung nodules.
```

Ground Truth Text: Two left lung nodules concerning for metastatic disease.
Left basilar opacity could represent atelectasis or consolidation.

Corrections:

```
{
  "0": {
    "corrections": "Two left lung nodules concerning for
```

FINERADSCORE

```

        metastatic disease.",
    "clinical severity": "Actionable nonurgent error",
    "comments": "Incorrect number of nodules, may affect longer
                term follow up",
    "error category": ["Incorrect location/position of finding",
                     "Incorrect severity of finding"]
},
"None": {
    "corrections": "Left basilar opacity could represent
                  atelectasis or consolidation.",
    "clinical severity": "Urgent error",
    "comments": "In the context of neutropenic fever, missed possible
                pneumonia is an urgent finding",
    "error category": ["Omission of finding"]
},
"1": {
    "corrections": "Not actionable",
    "clinical severity": "[delete]",
    "comments": "Repetitive",
    "error category": []
}
}
}

```

Example 3:

Generated Text:

```

[0] Stable position of endotracheal tube projects 2.2 cm above the carina.
[1] Minimal atelectasis at the right lung base. [2] Moderate cardiomegaly.
[3] Pulmonary edema.
[4] The presence of a minimal left pleural effusion cannot be excluded.

```

Ground Truth Text: Endotracheal tube projects approximately 2.2 cm above the carina. Minimal atelectasis at the left and right lung bases. Moderate cardiomegaly. The presence of a minimal left pleural effusion cannot be excluded.

Corrections:

```

{
    "0": {
        "corrections": "Endotracheal tube projects approximately 2.2 cm
                      above the carina.",
        "clinical severity": "Not actionable",
        "comments": "Report references a non existent comparison but this
                    does not affect management",
        "error category": ["Mention of comparison that is not present in

```

```

    the reference impression"]
  },
  "1": {
    "corrections": "Minimal atelectasis at the left and right lung
                    bases.",
    "clinical severity": "Not actionable",
    "comments": "Incorrect anatomy of atelectasis would not
                affect management",
    "error category": ["Incorrect location/position of finding"]
  },
  "3": {
    "corrections": "[delete]",
    "clinical severity": "Urgent error",
    "comments": "Incorrect call of pulmonary edema is a potentially
                urgent error if diuretics are started,
                for example",
    "error category": ["False prediction of finding"]
  }
}

```

Example 4:

Generated Text: [0] Endotracheal tube is in standard position.

[1] A nasogastric tube is
seen coursing into the stomach with tip in the stomach.
[2] Heart size is normal.
[3] Lungs are clear.
[4] No pleural effusion or pneumothorax. [5] Enteric tube
tip is in the stomach.

Ground Truth Text: Endotracheal tube is in standard position.

Heart size is normal.
Lungs are clear. No pleural effusion or pneumothorax.

Corrections:

```

{  "1": {
    "corrections": "[delete]",
    "clinical severity": "Actionable nonurgent error",
    "comments": "no gastric tube",
    "error category": ["False prediction of finding"]
  },
  "5": {
    "corrections": "[delete]",
    "clinical severity": "Actionable nonurgent error",
    "comments": "no enteric tube",
  }
}

```

FINERADSCORE

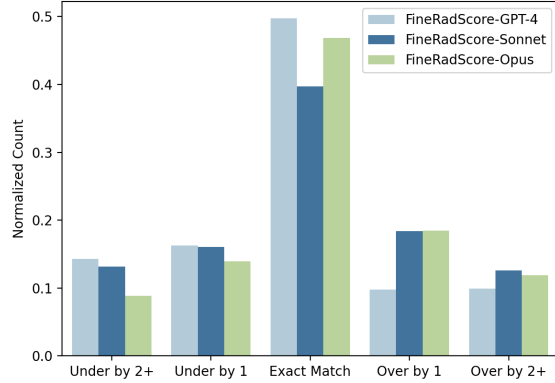


Figure 4: FineRadScore Clinical Severity Scores' Accuracies on ReFiSco-v0 with a Shuffled Ground Truth. FineRadScore's clinical severity scores line up closely with the ground truth radiologist's on ReFiSco-v0. FineRadScore-GPT-4 achieves the best performance in assigning clinical severity scores.

```
    "error category": ["False prediction of finding"]
  }
}
```

Example 5:

Generated Text: [0] The lungs are well expanded. [1] There is no pleural effusion or pneumothorax.
[2] The cardiomediastinal and hilar contours are unremarkable.

Ground Truth Text: The lungs are adequately inflated. The contours of the cardiomediastinal and hilar regions appear normal. There are no indications of pleural effusion or pneumothorax.

Corrections: {}

Appendix B. Clinical Severity Score Graphs

FINERADSCORE

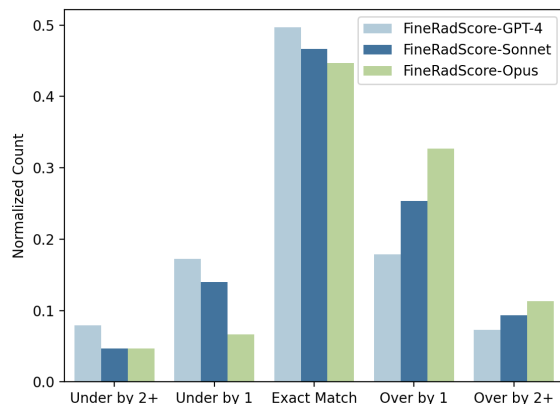


Figure 5: FineRadScore Clinical Severity Scores’ Accuracies on ReFiSco-v1 with a Paraphrased Ground Truth. FineRadScore’s clinical severity scores line up closely with the ground truth radiologist’s on ReFiSco-v1. FineRadScore-GPT-4 achieves the best performance in assigning clinical severity scores.

Appendix C. Radiologist Correction Type Graph

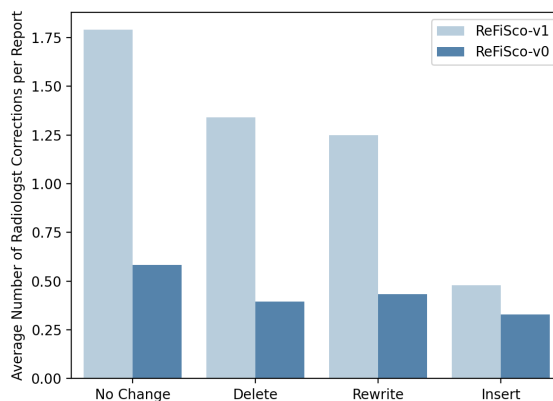


Figure 6: Radiologists leave most lines unchanged, and delete or rewrite lines at similar rates. Slightly under half of the reports contain at least one insertion by the radiologist.

Appendix D. FineRadScore Labeling ReXVal Error Categories

Overall, FineRadScore still falls short in labeling each correction it generates with the correct ReXVal category. To reiterate, the ReXVal error categories are as follows:

1. False prediction of finding

2. Omission of finding
3. Incorrect location/position of finding
4. Incorrect severity of finding
5. Mention of comparison that is not present in the reference impression
6. Omission of comparison describing a change from a previous study

For each report in the ReXVal dataset, for each of the categories above, we record whether or not the radiologist found at least one error in that category, as well as whether FineRadScore labeled at least one correction with that category. We use these numbers to compute precision and recall metrics across all six ReXVal categories. As shown in the table below, we see low precision numbers and high recall numbers. These results occur because FineRadScore is over-labeling corrections with these ReXVal error categories. It has a tendency to apply multiple labels to corrections, even when just one fits.

	Type 1	Type 2	Type 3	Type 4	Type 5	Type 6
Precision	32.62%	64.00%	33.33%	36.54%	40.91%	26.39%
Recall	100%	100%	100%	100%	100%	100%

Appendix E. FineRadScore Comment Quality

To evaluate the quality of FineRadScore-GPT-4 generated comments, we ask radiologists to evaluate if a subset of comments were 1) accurate and 2) helpful. The radiologists find 71.98% of the comments accurate and 79.74% helpful. We now offer a qualitative breakdown of the types of comments FineRadScore-GPT-4 produces. Some comments simply state that the line that needs to be corrected does not line up with the ground truth, which is technically true, but not informative. Some comments are correct, but referencing the wrong line in the report. Some comments are partially true so still helpful, which led to them getting a “not accurate” rating but still helpful.