

Maximizing NMR Sensitivity: A Guide to Receiver Gain Adjustment

Josh P. Peters^{1*}, Jan-Bernd Hövener¹, Andrey N. Pravdivtsev^{1*}

1. Section Biomedical Imaging, Molecular Imaging North Competence Center (MOIN CC), Department of Radiology and Neuroradiology, University Medical Center Kiel, Kiel University, Am Botanischen Garten 14, 24114, Kiel, Germany

*Corresponding Authors: josh.peters@rad.uni-kiel.de; andrey.pravdivtsev@rad.uni-kiel.de

Abstract

Nuclear magnetic resonance is revolutionizing in parallel with the development of methodology and all electronics and hardware around it. The primary objective of such modifications is to make the system more reliable and sensitive. Almost a decade ago, Bruker introduced an Avance NEO console with an extended dynamic range of the receiver gain (RG). The system offered a maximum sensitivity even at very low RG. This feature is handy in the case of hyperpolarization when too much signal can cause clipping of the free induction decay signal. We analyzed the RG profile for our spectrometer's ^1H , ^2H , ^{13}C , and ^{15}N nuclei and found non-monotonous sensitivity behavior characterized by signal-to-noise ratio (SNR). For example, ^{13}C SNR at RG = 18 is 5.5% higher than SNR at maximum RG = 101 and 51.6% higher than at RG = 20.2. Having found this, we advertise here that other NMR users using any spectrometer should calibrate their SNR(RG) and find the RG value right for their experiments, as automatic RG adjustment does not know about the actual performance of each system. We offer a tool and strategy for such an optimization and exemplify the benefits of using SNR-adjusted RG value with hyperpolarized [$1\text{-}^{13}\text{C}$]pyruvate measurements.

Introduction

Nuclear magnetic resonance is a universal analytical method¹. The development of stronger magnets, less noisy electronics, and more sensitive probes, including cryoprobes, further extends the scope of the technique^{2,3}. Still, the sensitivity of NMR is a hurdle that incentivizes members of the hyperpolarization community to develop various signal enhancement methods⁴⁻⁶.

Maximizing the signal-to-noise ratio (SNR) is a nontrivial task. Typically, it is advised to use automatic receiver gain adjustment, which maximizes RG while avoiding the signal overflow and clipping of the free induction decay (FID, **Figure 1**). The clipping leads to spectral artifacts, which would make the spectrum useless. However, in the case of hyperpolarization, automated RG adjustment is typically not possible in the first place, as the signal is enhanced only transiently and rapidly decays. Therefore, sufficiently low RG and

excitation pulses are used to avoid the overflow. However, this approach may be too pessimistic and uninformed, resulting in a loss of SNR for the sake of clipping-free FID.

For a long time, Bruker introduced a dynamic receiver gain that provides sufficiently high sensitivity even at low RG values available from Avance NEO console. Having such a console and studying real-time metabolomics using the dissolution dynamic nuclear polarization method for signal enhancement (dDNP)⁵, we wanted to find the optimal RG and excitation angle value apriori, knowing the polarization and concentration of polarization of the system under investigation.

To our surprise, SNR was not monotonous as a function of RG (**Figure 2-SNR**), indicating that standard automated RG adjustment routine can not adequately set this value even for standard experiments without hyperpolarization. Below, we will discuss how to find SNR figures as a function of

RG (i), optimal RG, and excitation angle parameters for thermal- and hyperpolarization experiments (ii) and exemplify the advantages of calibrating the RG by recording hyperpolarized $[1-^{13}\text{C}]$ pyruvate (iii). This calibration procedure does not require much time and offers free SNR gain.

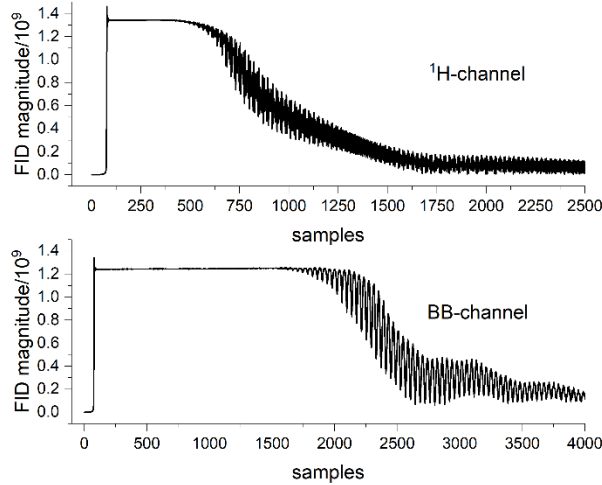


Figure 1: FID magnitude of pure H₂O (top) and 99 mM hyperpolarized $[1-^{13}\text{C}]$ pyruvate (bottom) measured such that the receiver overflows. The natural shape of the FID is altered by clipping the signal at the beginning – flat signal at maximum ADC signal. The maximum amplitude of the ADC for each channel can be obtained from this ADC overflow experiment. The respective values are 1.344×10^9 (^1H -channel) and 1.244×10^9 (BB-channel).

Materials and Methods

NMR system: All experiments were performed on a 9.4-T wide bore NMR (WB400, Avance NEO, Bruker) with a 5-mm narrow bore BBFO probe. Standard 5 mm economy NMR tubes filled to 550 μL were used as sample tubes.

^1H SNR calibration sample: 9.93 mg H₂O, 580.87 mg D₂O (00506, Deutero GmbH), and 13.34 mg (10 μL) Gadolinium (Gd, gadobutrol Gd-DO3A-butrol, 1 mmol/ml; Gadovist, Bayer). 10 μL Gd contain about 7.29 mg of H₂O, in total this converts to 17.17 μL H₂O and 525.91 μL D₂O, leading to a 1.728 M H₂O and a 3.456 M ^1H concentration.

^2H SNR calibration sample: 527.28 mg H₂O, 11.09 mg D₂O, and 13.64 mg (10 μL) Gd. This converts to 528.87 μL H₂O, 10.04 μL D₂O and

10 μL Gd. This leads to a 1.009 M D₂O and 2.018 M ^2H concentration.

^{13}C SNR calibration sample: 400 μL H₂O, 100 μL D₂O, and 50 μL $[1-^{13}\text{C}]$ pyruvate (677175, Sigma-Aldrich). This leads to a 1.308 M $[1-^{13}\text{C}]$ pyruvate and 1.308 M ^{13}C concentration.

^{15}N SNR calibration sample: 494 mg H₂O, 55.07 mg D₂O, 55.08 mg $^{15}\text{NH}_4$ (299251, Sigma-Aldrich), and 10 μL Gd. This converts to 495.5 μL H₂O and 49.87 μL D₂O. This leads to a 2.059 M $^{15}\text{NH}_4$ and 2.059 M ^{15}N concentration.

Measurement of SNR as a function of RG: The NMR signal amplitude is affected by several parameters:

$$\text{signal} = \text{signal}(RG, \dots) = A \cdot RG \cdot \sin(\alpha) \cdot P \cdot C \quad \text{eq 1}$$

where RG is the value of linear receiver gain, α is the angle of the excitation pulse, P is the polarization value of nuclear spin of interest, C is the spin concentration, and A is a hardware coefficient independent of RG , α , P and C which is easy to find experimentally.

The ADC system has a threshold for the maximum signal intensity S_m that puts a restriction on the allowed signal:

$$\text{signal} \leq S_m \quad \text{eq 2}$$

The noise level is a nontrivial function of many parameters, including RG value:

$$\text{noise} = \text{noise}(RG, \dots) \quad \text{eq 3}$$

The signal-to-noise ratio, SNR, then can be estimated as

$$\text{SNR}(RG) = \frac{\text{signal}(RG, \dots)}{\text{noise}(RG, \dots)} = \frac{A \cdot RG \cdot \sin(\alpha) \cdot P \cdot C}{\text{noise}(RG, \dots)} \quad \text{eq 4}$$

Since some parameters of eq 4 are unknown (noise, A) it is necessary to calculate the reference $\text{SNR}_{\text{ref}}(RG)$ with all the known parameters such that this can be used to estimate the SNR of other nuclear spin resonances:

$$\text{SNR}(RG) = \text{SNR}_{\text{ref}}(RG) \cdot \frac{\sin(\alpha)}{\sin(\alpha_{\text{ref}})} \cdot \frac{P}{P_{\text{ref}}} \cdot \frac{C}{C_{\text{ref}}} \quad \text{eq 5}$$

We used one of the calibration samples described above and the following protocol to measure $\text{SNR}_{\text{ref}}(RG)$:

1. Place the NMR tube with 550 μL of calibration sample into the NMR probe and wait for temperature equilibration.
2. Adjust probe resonance and shims and calibrate for 90° RF pulse.
3. Use the standard excitation and measurement protocol α -FID to measure FIDs as a function of RG for the known other parameters: α , P , C . Ensure that the time between experiments is long enough, at least $3 \times T_1$ (consider adding Gd). Repetition time TR was 30 s for ^1H , 20 s for ^2H , 300 s for ^{13}C , and 80 s for ^{15}N (e.g., $T_1^{15\text{N}} = 2.1$ s). We used 2 averaging scans for ^{13}C , 8 for ^1H and ^2H , and 16 for ^{15}N , and 1 dummy scan for ^{13}C or 4 for ^1H , ^2H , and ^{15}N .

Having obtained the FIDs, we used MestReNova software (version 14.2.0-26256, 2020) to evaluate the signal intensities and SNR. SNR was calculated using the standard tool “*SNR Graph*” where you define the signal region, and the tool automatically finds the noise area and calculates the SNR of the system. The tool calculates the SNR using a ratio between those parameters: $\text{SNR} = Y / \text{noise}$, where Y is the amplitude of the peak and noise is the root-mean-square deviation calculated in a signal-free area of the spectrum. Using this protocol, we measured the reference signal and SNR (**Figure 2**).

Calculation of optimal RG: It is a classical constraint-satisfaction problem to find the maximum for SNR (eq. 5) as a function of many parameters with constraints. In the experiment, essentially, only two parameters that are nonlinearly dependent on SNR can change: RG and α . The rest are predetermined or can be treated then as semifixed.

In this case, the complete optimization task can be written as:

$$\begin{cases} \text{SNR}(RG) \rightarrow \max \\ \text{signal} \leq S_m \\ 0 < \alpha \leq \alpha_m \leq 90^\circ \end{cases} \quad \text{eq 6}$$

This ensures that the signal is below the threshold signal value, S_m , so ADC does not overflow (**Figure 1**, condition 2). In addition, one can put supplementary constraints on the excitation angle (condition 3): this is especially important in hyperpolarization experiments where signal is measured using small α . Hence, using this eq 6, one can maximize SNR for given P , C , and maximum allowed α_m for the experiment.

Because there are only two variables, α and RG , which one can typically change, and the experiment defines the P and C , one can quickly maximize SNR under the constraints of eq 6 numerically. We propose the following approach, which is visually easy to understand. The algorithm is realized for our system using Python script and is available as a Supporting Material.

The approach consists of the following five steps:

1. Experimentally measure $\text{SNR}_{\text{ref}}(RG)$ for constant α_{ref} , P_{ref} , C_{ref} . We used thermally polarized reference samples for that, as described before.
2. Specify required S_m and α_m .
3. Calculate $\text{SNR}(RG, \alpha)$ and $\text{signal}(RG, \alpha)$ for available RG region (for NEO it is from 0.25 to 101) and $0 \leq \alpha \leq \alpha_m$.
4. Find the mask when $\text{signal}(RG, \alpha) \leq S_m$.
5. Apply the mask for $\text{SNR}(RG, \alpha)$ and find the maximum.

Alternatively to the numerical approach, one can use the Karush–Kuhn–Tucker method for finding the function's maximum with parameter constraints.

dDNP experiment:

Following established protocol⁷, the pyruvate stock sample was prepared by mixing about 25 mg of trityl radical (AH111501, Polarize) and $[1\text{-}^{13}\text{C}]$ pyruvate (677175, Sigma-Aldrich). This resulted in 31 mM trityl and 14 M pyruvate concentrations in the stock sample, and 20 μL were used for hyperpolarization. After dissolution with 3.7 mL of superheated ($\sim 200^\circ\text{C}$, 11 bar) dissolution medium, the sample was transferred to an NMR system and detected after about 18 s. Using such polarized pyruvate, we carried out two experiments.

First, we measured the threshold for the signal on the BB channel. The hyperpolarized sample was placed in the NMR and the signal was clipped (see below). The polarization at the measuring site was about 38%.

Second, we measured the NMR signal every 7-10 s using $\alpha = 5^\circ$ and alternating RG between 0.25 and 18. The repetition time between experiments was not constant because of the internal delay of the console needed to compile and execute individual experiments. The polarization at the measuring site was 32.9%.

Method to find S_m :

To find the signal threshold S_m we measured NMR signals of pure H_2O (1H TxRx channel) with $\alpha = 90^\circ$ or hyperpolarized $[1-^{13}C]$ pyruvate (BB TxRx channel) with $\alpha = 40^\circ$ and the maximum $RG = 101$ and observed the maximum signal in the FID, indicating the threshold of the spectrometer's ADC (**Figure 1**). The maximum value for 1H was $S_m \sim 1.34 \cdot 10^9$ and for ^{13}C it was $S_m \sim 1.24 \cdot 10^9$.

Results.

SNR as a function of RG

To find the actual performance of the NMR system, we acquired NMR spectra for four reference samples as a function of RG (**Figure 2**). Each sample resulted in a dominant, intense NMR resonance.

For the BB channel, we found about a 40% dip in SNR figures at RG between 20 and 60 compared to $RG = 18$, although the signal increased linearly throughout the range.

The highest and somewhat similar SNR for 2H , ^{13}C , and ^{15}N was measured at $RG = 8$ to 18 with deviations in SNR of only 1-2%. The SNR for 1H was still slowly rising, reaching its maximum at $RG = 101$. This finding shows that the most sensitive RG for X nuclei is about 18 and should be used in most cases without any automatic RG adjustment. The lower RG values are needed when the signal exceeds the threshold. For the 1H , however, a small gain in SNR is still present above $RG = 18$.

This finding revealed that maximizing RG is not the right strategy for improving SNR.

To check for the linearity between RG and signal, we used a formula eq 7, which shows the deviation from the $signal_{max} = signal(RG = RG_{max} = 101)$ and its respective extrapolation down to $RG = 0.25$:

$$signal_{deviation}(RG) = \frac{signal(RG)}{RG} \cdot \frac{RG_{max}}{signal_{max}} - 1 \quad \text{eq. 7}$$

This revealed that the 1H and 2H measurements are quite linear (maximum deviation 1-2%), while ^{13}C and ^{15}N have quite substantial deviations, especially for low RG values (about 8% and 30%, respectively, **Figure 2**-third column). Despite the signal and SNR being about 100 times lower

compared to 1H and 2H , still, the SNR for ^{13}C and ^{15}N is about 50 even for the lowest RG , indicating that observations can not be explained by low sensitivity. However, it may suggest that lower signals are more prone to experience RG non-linearity which should be checked experimentally.

Setting optimal excitation for your hyperpolarization sample

Having found the $SNR(RG)$ function (**Figure 2**) and measuring the FID clipping threshold for each nucleus: $S_m \sim 1.34 \cdot 10^9$ for 1H and $S_m \sim 1.24 \cdot 10^9$ for BB channels (**Figure 1**), one can estimate the signal intensity and SNR for a known sample volume and composition for any parameters: RG , α , C , and P .

For 1H and 2H , we checked when pure water (H_2O or D_2O), the most typical solvent, can be measured without ADC overflow. In both cases, thermal polarization was used, and concentrations $[H_2O]$ (**Figure 3- 1H**) and $[D_2O]$ (**Figure 3- 2H**) were set to 55.51 M and 55.21 M, respectively. The corresponding number of spins is double the concentration. We observe that in both cases, the ADC can be easily clipped, so α and RG should be chosen with caution.

Note that when one aims to maximize SNR without restraints for the excitation angle, maximizing α is more beneficial compared to maximizing RG (**Figure S1**, SI). For water, the maximum SNR was obtained with $\alpha = 90^\circ$ and RG of 4.9 and 12.8, respectively.

In hyperpolarization experiments, the C and P values are typically fixed, and only α and RG are free parameters. We calculated the signal and SNR maps as a function of RG and α for a typical ^{13}C sample: concentrations of single ^{13}C labeled pyruvate or fumarate of 90 mM with a polarization of 35% (**Figure 3- ^{13}C**), and ^{15}N sample 40 mM and a polarization of 15% - similar to values achieved with $[1-^{15}N]$ nicotinamide previously⁸ (**Figure 3- ^{15}N**). The flipping angle α in hyperpolarization metabolic experiments is typically low ($\sin(\alpha) \sim \alpha$) to be able to track the polarization over a long period. Here (**Figure 3- ^{13}C** , **- ^{15}N**), they have been limited to 5° (^{13}C) and 10° (^{15}N).

For the given conditions, we calculated the maximum product of $\sin(\alpha) \cdot RG$ that can be used without clipping, resulting in values 4.99 for 1H ,

24.82 for ^2H , 4.14 for ^{13}C , and 107.15 for ^{15}N . This means that, e.g. for hyperpolarized ^{13}C with a 90° pulse, only $RG \leq 4.14$ is possible to avoid clipping, while for ^{15}N , any combination of RG (from 0.25 to 101) and flipping angle (from 0 to 90°) can be used.

Considering also the nonlinear $SNR(RG)$ dependence (**Figure 2**), the maximum SNR is obtained for all X-nuclei and $\alpha < 5^\circ$ at $RG =$

8 to 18, where $SNR(RG)$ dependence reached its maximum before going down.

Note that here, we use an exact condition for the signal overflow: $\text{signal} \geq S_m$, however, it is maybe practical to use a more conservative condition that leaves some room for errors without actual signal overflow. E.g., signal must not exceed half of the threshold value: $\text{signal} \geq 0.5 \cdot S_m$.

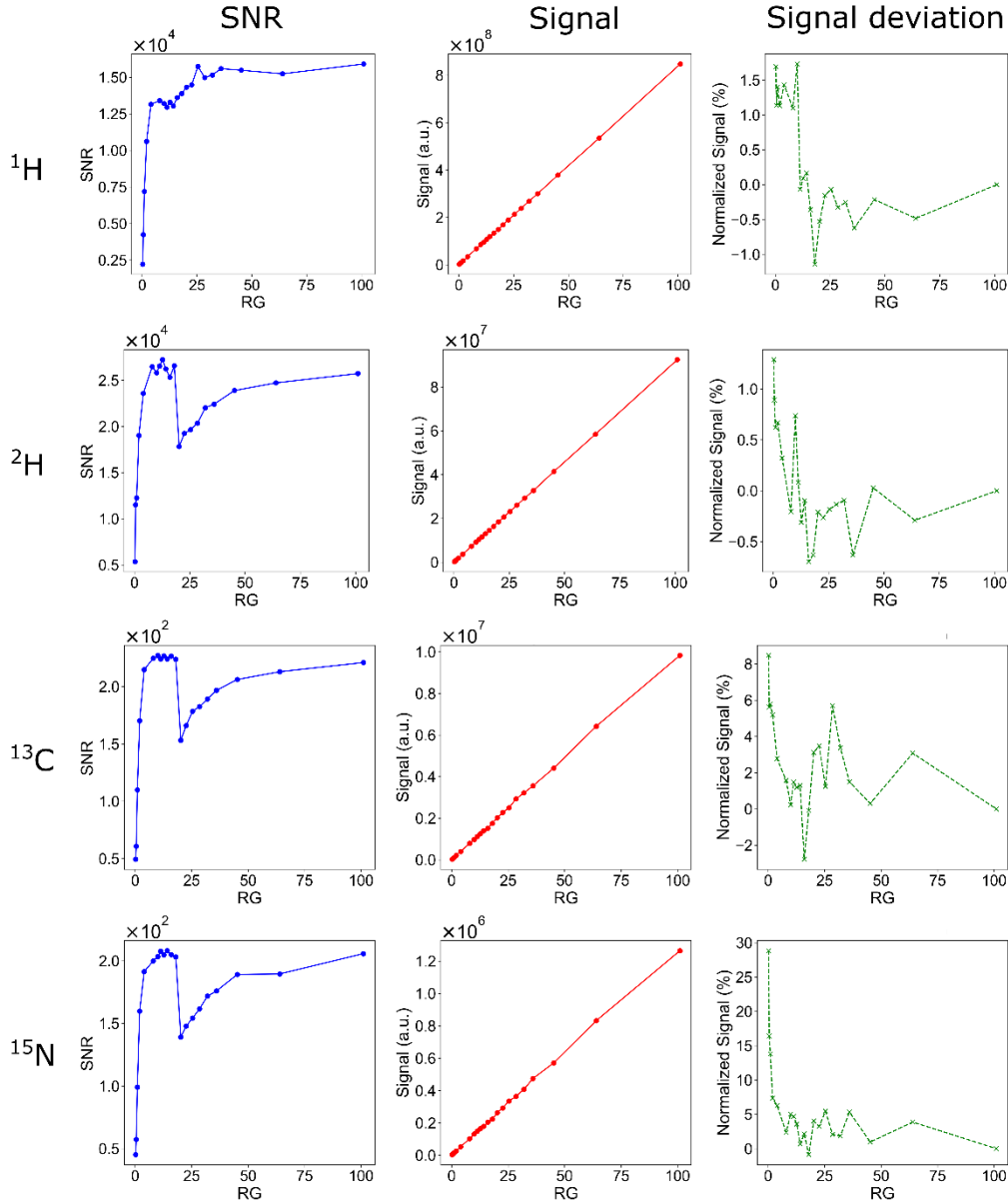


Figure 2: SNR (first column), corresponding signals (middle column), and deviation of signal amplitude from the linear signal gain (right column) for ^1H , ^2H , ^{13}C , and ^{15}N reference samples as a function of RG . The reference sample composition is given in the methods section. Note the maximum and a dip in $SNR(RG)$ after $RG = 18$ and the prolonged rate of SNR recovery for all samples except ^1H .

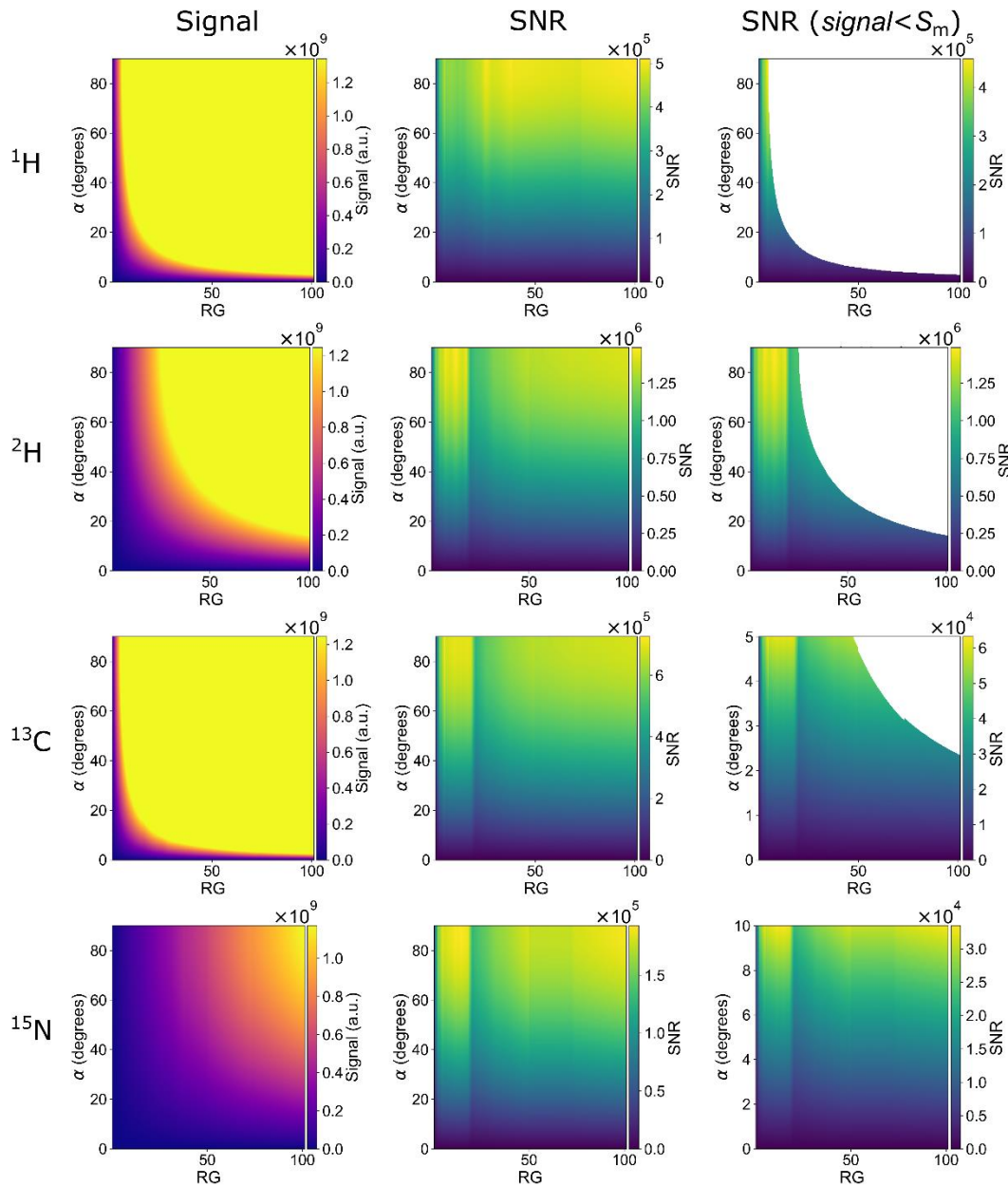


Figure 3: Signal map (left column), SNR map (middle column), and clip-free SNR map (right column) calculation example. The following parameters were used in simulations, which should reflect possible experimental conditions: (^1H) a pure H_2O sample at thermal polarization was assumed with $[\text{H}_2\text{O}] = 55.51 \text{ M}$, (^2H) a pure D_2O sample with $[\text{D}_2\text{O}] = 55.21 \text{ M}$, (^{13}C) a 90 mM ^{13}C with 35% polarization, and (^{15}N) a 40 mM ^{15}N with 15% polarization. The maps were calculated by extrapolating signal(RG) and SNR(RG) reference maps (**Figure 2**) for missing RG values between 0.25 and 101. Values when signal $> S_m$ were removed from the SNR maps (third column, white areas), leaving behind only allowed (signal $< S_m$, when ADC is not overflowed) combinations of α and RG.

Experimental evidence of benefits of optimizing RG in hyperpolarization experiment

We measured the NMR signal of hyperpolarized to 33% $[\text{1-}^{13}\text{C}]\text{pyruvate}$ every 7-10 s using $\alpha = 5^\circ$ and

alternating RG between 0.25 and 18 (**Figure 4**). Using this approach we were able to compare the sensitivity gain in the same hyperpolarization experiment. The difference in SNR by a factor of 5

is apparent, clearly showcasing the gained sensitivity. The signal was normalized using the ^{13}C calibration experiment (**Figure 2**), where a signal ratio of 66.33 (instead of $18 / 0.25 = 72$ expected for perfectly linear RG) was observed between signals measured at RG 0.25 and 18. Following normalization, the monoexponential decay inherent to hyperpolarized signals is apparent despite the alternating RGs.

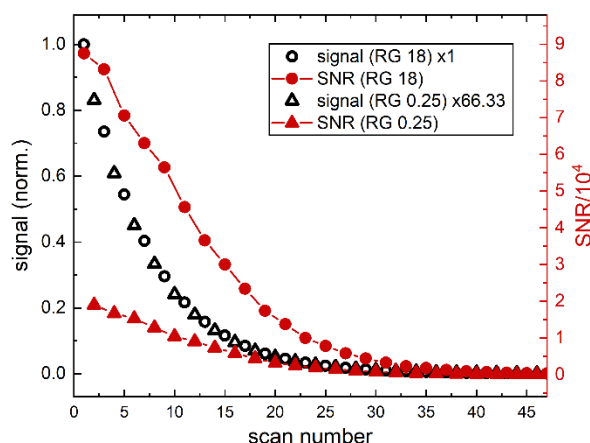


Figure 4: Signal and SNR of hyperpolarized $[1-^{13}\text{C}]$ pyruvate as a function of scan number acquired using flipping angle, $\alpha = 5^\circ$, and alternating RG between 0.25 and 18. The signal of RG 0.25 was multiplied by calibrated factor 66.33 (instead of linear 72, **Figure 2) to bring it to the same level as the signal obtained with RG 18. The monotonous decay of the signal and SNR is clearly visible. The repetition time is about 7 to 10 s. The difference in SNR by a factor of about 5 between the two RG settings is apparent, which is similar to the ratio in SNR obtained from the calibration experiment equal to ~ 4.8 .**

Conclusion

In one of our previous studies⁹ we used RG of 10 for measuring hyperpolarized ^{13}C signal. Here, we found that SNR at RG of 10 is 1.5% lower than at the most sensitive RG of 18 and 2.8% higher than at maximal RG of 101. Our measurements here reveal that we were still below the FID clipping threshold and very close to optimal RG in those experiments. By more careful adjustment of SNR, using 18 instead of 10 or 101, we could have improved SNR by 1.5% and 1.3%, respectively, meaning that on our system, there are specific RG ranges where SNR is very close to the maximum

despite different signal amplitudes. At the same time, however, if we chose an RG of 20.2, the SNR would be about half of the SNR gained at RG 10, 18, or 101. Notably, improving SNR by only 50% gives a factor of 2 acceleration when signal averaging is required.

The $\text{SNR}(\text{RG})$ profile reported here was measured only for one available system. However, the central message of the manuscript does not depend on the specific system one uses: Although it is known that the modern NMR system is very reliable, it is still worth checking your $\text{SNR}(\text{RG})$ profile and setting the parameters right for your system to gain maximum SNR in hyperpolarization or conventional NMR experiments.

Acknowledgements

We acknowledge funding from German Federal Ministry of Education and Research (BMBF) within the framework of the e:Med research and funding concept (01ZX1915C, 03WIR6208A hyperquant), DFG (PR 1868/3-1, PR 1868/5-1, HO-4602/2-2, HO-4602/3, HO-4602/4, EXC2167, FOR5042, TRR287). MOIN CC was founded by a grant from the European Regional Development Fund (ERDF) and the Zukunftsprogramm Wirtschaft of Schleswig-Holstein (Project no. 122-09-053).

Supporting Materials

FIDs, and additional spectra are available in supporting materials (.pdf). Raw data are available in Zenodo repository: DOI: 10.5281/zenodo.13795731.

Author Information

Josh Peters, ORCID 0000-0003-1019-4067
Prof. Jan-Bernd Hövener: ORCID 0000-0001-7255-7252;
Dr. Andrey N. Pravdivtsev, ORCID 0000-0002-8763-617X

Author Contributions

ANP, JBH: conceptualization; ANP, JP: investigation, analysis, writing – original draft; ANP and JBH: supervision, funding acquisition. All authors contributed to discussions and

interpretation of the results and have approved the final version of the manuscript.

References

1. Sapir G, Steinberg DJ, Aqeilan RI, Katz-Brull R. Real-Time Non-Invasive and Direct Determination of Lactate Dehydrogenase Activity in Cerebral Organoids—A New Method to Characterize the Metabolism of Brain Organoids? *Pharmaceuticals*. 2021;14(9):878. doi:10.3390/ph14090878
2. Sharma R, Sharma A. 21.1 Tesla Magnetic Resonance Imaging Apparatus and Image Interpretation: First Report of a Scientific Advancement. *Recent Pat Med Imaging*. 2011;1(2):89-105.
3. Kovacs H, Moskau D, Spraul M. Cryogenically cooled probes—a leap in NMR technology. *Prog Nucl Magn Reson Spectrosc*. 2005;46(2):131-155. doi:10.1016/j.pnmrs.2005.03.001
4. Hövener J, Pravdivtsev AN, Kidd B, et al. Parahydrogen-Based Hyperpolarization for Biomedicine. *Angew Chem Int Ed*. 2018;57(35):11140-11162. doi:10.1002/anie.201711842
5. Ardenkjaer-Larsen JH, Fridlund B, Gram A, et al. Increase in signal-to-noise ratio of > 10,000 times in liquid-state NMR. *Proc Natl Acad Sci*. 2003;100(18):10158-10163. doi:10.1073/pnas.1733835100
6. Eills J, Budker D, Cavagnero S, et al. Spin Hyperpolarization in Modern Magnetic Resonance. *Chem Rev*. 2023;123(4):1417-1551. doi:10.1021/acs.chemrev.2c00534
7. Ferrari A, Peters J, Anikeeva M, et al. Performance and reproducibility of ¹³C and ¹⁵N hyperpolarization using a cryogen-free DNP polarizer. *Sci Rep*. 2022;12(1):11694. doi:10.1038/s41598-022-15380-7
8. Peters JP, Brahms A, Janicaud V, et al. Nitrogen-15 dynamic nuclear polarization of nicotinamide derivatives in biocompatible solutions. *Sci Adv*. 2023;9(34):eadd3643. doi:10.1126/sciadv.add3643
9. Peters JP, Assaf C, Haj Mohamad F, et al. Yeast solutions and hyperpolarization enable real-time observation of metabolized substrates even at natural abundance. *Unpublished*.

Supporting information for:

Maximizing NMR Sensitivity: A Guide to Receiver Gain Adjustment

Josh P. Peters^{1*}, Jan-Bernd Hövener¹, Andrey N. Pravdivtsev^{1*}

1. Section Biomedical Imaging, Molecular Imaging North Competence Center (MOIN CC), Department of Radiology and Neuroradiology, University Medical Center Kiel, Kiel University, Am Botanischen Garten 14, 24114, Kiel, Germany

*Corresponding Authors: josh.peters@rad.uni-kiel.de; andrey.pravdivtsev@rad.uni-kiel.de

Contents

Constant signal experiment	2
¹ H data	3
² H data	5
¹³ C data	7
¹⁵ N data	9

Constant signal experiment

Constant signal experiment as a function of RG: We did one experiment where we kept a constant $\text{signal}(\text{RG}) = \text{const.}$ According to eq 1, this can be achieved by setting $\text{RG} \cdot \sin(\alpha)$ to a constant value. If one sets $\alpha = 90^\circ$ for the lowest RG of 0.25, then

$$\sin(\alpha) = \frac{0.25}{\text{RG}} \quad \text{eq S1}$$

gives the α angle as a function of RG, providing a constant signal. Using this approach, we acquired an RG sweep with varied flip angles from 90° (RG 0.25) to 0.142° (RG 101) to keep a constant signal amplitude. The highest SNR was obtained for the combination of RG 0.25 and 90° FA and subsequently fell monotonously. This indicates that for a given signal threshold, it is best first to maximize the flip angle within the constraints of the experiment and then maximize RG as a second measure. However, maximizing the flip angle is typically not feasible, especially during a hyperpolarization experiment. In this case, the optimization procedure from the main text applies.

Note that the signal decreases from $0.85 \cdot 10^5$ a.u. to $0.60 \cdot 10^5$ a.u. between RG 0.25 to 45.2, which is consistent with the finding from the RG sweeps in **Figure 2** in the main text where the signal at lower RG values is over-proportionally high.

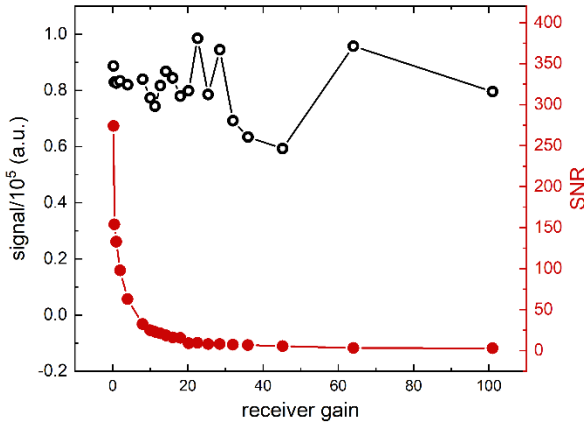


Figure S1: Constant signal experiment on the ^{15}N -sample showcases the linearity of the signal and calibration of the 90° pulse. The signal gets noisier for higher RGs, because the flipping angle gets as low as 0.142° at RG 101. It also demonstrates that for a given threshold (here $\text{signal} \approx 0.85 \cdot 10^5$ a.u.) one should maximise first the flipping angle and then RG if maximum SNR is the goal. It is apparent, similar to Figure 3- ^{15}N from the main text, that the signal is higher for lower RG values if $\sin(\alpha) \cdot \text{RG}$ is kept constant. In theory, the signal should stay constant if the RG is linear, however here, it decreases from $0.85 \cdot 10^5$ a.u. (RG 0.25) to $0.60 \cdot 10^5$ a.u. (RG 45.2), a decrease by 30%.

¹H data

Table S1: Signal-to-noise-ratio (SNR) and maximum FID magnitude values measured with ¹H sample depending on receiver gain (RG).

RG	SNR	max FID signal (a.u.)
0.25	2.21E+03	2.13E+06
0.5	4.23E+03	4.24E+06
1	7.19E+03	8.51E+06
2	1.06E+04	1.70E+07
4	1.32E+04	3.40E+07
8	1.34E+04	6.79E+07
10	1.32E+04	8.54E+07
11.3	1.30E+04	9.48E+07
12.7	1.33E+04	1.07E+08
14.2	1.30E+04	1.19E+08
16	1.36E+04	1.34E+08
18	1.39E+04	1.49E+08
20.2	1.43E+04	1.69E+08
22.6	1.45E+04	1.89E+08
25.4	1.57E+04	2.13E+08
28.5	1.50E+04	2.38E+08
32	1.52E+04	2.68E+08
36	1.56E+04	3.00E+08
45.2	1.55E+04	3.78E+08
64	1.52E+04	5.34E+08
101	1.59E+04	8.47E+08

Data generated from analysis: /1H

+++++

Analysis settings:

Alpha reference: 90°

Averages reference: 8

Concentration reference: 3.456 M

Polarization reference: 0.0030981%

Concentration experiment: 111.02 M

Polarization experiment: 0.0030981%

Maximum ADC amplitude: 1.344E+09 a.u.

Extracted RG array length: 21

Extracted SNR array length: 21

Extracted signal array length: 21

Lowest RG value which yields 95% of the maximum SNR: 32

Maximum SNR in the masked plot: 459230

Corresponding RG values of max SNR: 4.89379

Corresponding flipping-angle (degrees) of max SNR: 90

Maximum product of sin(alpha) and RG with chosen concentration (111.02 M) and polarization (0.0030981%) before clipping: 4.98604

Your molar concentration in use was: 0.00343951

Maximum flip angle for RG 32 and molar concentration 0.00343951: 8.83768°

Maximum valid molar concentration for RG 32 and 90°: 0.00053728

+++++

^2H data

Table S2: Signal-to-noise-ratio (SNR) and maximum FID magnitude values measured with ^2H sample depending on receiver gain (RG):

RG	SNR	max FID signal (a.u.)
0.25	5.34E+03	2.32E+05
0.5	1.15E+04	4.62E+05
1	1.23E+04	9.22E+05
2	1.90E+04	1.84E+06
4	2.36E+04	3.68E+06
8	2.65E+04	7.31E+06
10	2.58E+04	9.23E+06
11.3	2.65E+04	1.04E+07
12.7	2.72E+04	1.16E+07
14.2	2.62E+04	1.30E+07
16	2.53E+04	1.46E+07
18	2.66E+04	1.64E+07
20.2	1.78E+04	1.85E+07
22.6	1.92E+04	2.06E+07
25.4	1.96E+04	2.32E+07
28.5	2.04E+04	2.61E+07
32	2.20E+04	2.93E+07
36	2.24E+04	3.28E+07
45.2	2.39E+04	4.14E+07
64	2.47E+04	5.85E+07
101	2.57E+04	9.25E+07

Data generated from analysis: /2H

+++++

Analysis settings:

Alpha reference: 90°

Averages reference: 8

Concentration reference: 2.018 M

Polarization reference: 0.000475899%

Concentration experiment: 110.42 M

Polarization experiment: 0.000475899%

Maximum ADC amplitude: 1.244E+09 a.u.

Extracted RG array length: 21

Extracted SNR array length: 21

Extracted signal array length: 21

Lowest RG value which yields 95% of the maximum SNR: 14.2

Maximum SNR in the masked plot: 1.48732E+06

Corresponding RG values of max SNR: 12.768

Corresponding flipping-angle (degrees) of max SNR: 90

Maximum product of sin(alpha) and RG with chosen concentration (110.42 M) and polarization (0.000475899%) before clipping: 24.8207

Your molar concentration in use was: 0.000525487

Maximum flip angle for RG 14.2 and molar concentration 0.000525487: 90°

Maximum valid molar concentration for RG 14.2 and 90°: 0.000919346

+++++

¹³C data

Table S3: Signal-to-noise-ratio (SNR) and maximum FID magnitude values measured with ¹³C sample depending on receiver gain (RG).

RG	SNR	max FID signal (a.u.)
0.25	4.94E+01	2.64E+04
0.5	6.07E+01	5.14E+04
1	1.10E+02	1.03E+05
2	1.70E+02	2.05E+05
4	2.15E+02	4.00E+05
8	2.25E+02	7.91E+05
10	2.27E+02	9.76E+05
11.3	2.24E+02	1.12E+06
12.7	2.27E+02	1.25E+06
14.2	2.24E+02	1.40E+06
16	2.27E+02	1.51E+06
18	2.24E+02	1.75E+06
20.2	1.53E+02	2.03E+06
22.6	1.66E+02	2.28E+06
25.4	1.79E+02	2.50E+06
28.5	1.83E+02	2.93E+06
32	1.89E+02	3.22E+06
36	1.97E+02	3.56E+06
45.2	2.06E+02	4.41E+06
64	2.13E+02	6.42E+06
101	2.21E+02	9.83E+06

Data generated from analysis: /13C

+++++

Analysis settings:

Alpha reference: 90°

Averages reference: 2

Concentration reference: 1.308 M

Polarization reference: 0.00077917%

Concentration experiment: 0.09 M

Polarization experiment: 35%

Maximum ADC amplitude: 1.244E+09 a.u.

Extracted RG array length: 21

Extracted SNR array length: 21

Extracted signal array length: 21

Lowest RG value which yields 95% of the maximum SNR: 101

Maximum SNR in the masked plot: 63251.6

Corresponding RG values of max SNR: 4.89379

Corresponding flipping-angle (degrees) of max SNR: 5

Maximum product of $\sin(\alpha)$ and RG with chosen concentration (0.09 M) and polarization (35%) before clipping: 4.1353

Your molar concentration in use was: 0.0315

Maximum flip angle for RG 101 and molar concentration 0.0315: 2.34469°

Maximum valid molar concentration for RG 101 and 90°: 0.0147979

+++++

¹⁵N data

Table S4: Signal-to-noise-ratio (SNR) and maximum FID magnitude values measured with ¹⁵N sample depending on receiver gain (RG).

RG	SNR	max FID signal (a.u.)
0.25	4.52E+01	4.03E+03
0.5	5.73E+01	7.29E+03
1	9.91E+01	1.42E+04
2	1.60E+02	2.69E+04
4	1.91E+02	5.32E+04
8	2.00E+02	1.02E+05
10	2.03E+02	1.31E+05
11.3	2.08E+02	1.48E+05
12.7	2.05E+02	1.65E+05
14.2	2.08E+02	1.79E+05
16	2.05E+02	2.05E+05
18	2.03E+02	2.23E+05
20.2	1.39E+02	2.63E+05
22.6	1.48E+02	2.92E+05
25.4	1.54E+02	3.35E+05
28.5	1.62E+02	3.64E+05
32	1.72E+02	4.08E+05
36	1.76E+02	4.75E+05
45.2	1.89E+02	5.71E+05
64	1.90E+02	8.32E+05
101	2.06E+02	1.26E+06

Data generated from analysis: /15N

+++++

Analysis settings:

Alpha reference: 90°

Averages reference: 16

Concentration reference: 2.059 M

Polarization reference: 0.00031414%

Concentration experiment: 0.04 M

Polarization experiment: 15%

Maximum ADC amplitude: 1.244E+09 a.u.

Extracted RG array length: 21

Extracted SNR array length: 21

Extracted signal array length: 21

Lowest RG value which yields 95% of the maximum SNR: 8

Maximum SNR in the masked plot: 33507.7

Corresponding RG values of max SNR: 14.1814

Corresponding flipping-angle (degrees) of max SNR: 10

Maximum product of $\sin(\alpha)$ and RG with chosen concentration (0.04 M) and polarization (15%) before clipping: 107.151

Your molar concentration in use was: 0.006

Maximum flip angle for RG 8 and molar concentration 0.006: 90°

Maximum valid molar concentration for RG 8 and 90°: 0.452103

+++++