

From Bedside to Desktop: A Data Protocol for Normative Intracranial EEG and Abnormality Mapping

Heather Woodhouse¹, Sarah J. Gascoigne¹, Gerard Hall¹, Callum Simpson¹, Nathan Evans¹, Gabrielle M. Schroeder¹, Peter N. Taylor^{1,2,3} and Yujiang Wang^{1,2,3,*}

¹CNNP Lab (www.cnnp-lab.com), Interdisciplinary Computing and Complex BioSystems Group, School of Computing, Newcastle University, Newcastle upon Tyne, United Kingdom

²Faculty of Medical Sciences, Newcastle University, Newcastle upon Tyne, United Kingdom

³UCL Queen Square Institute of Neurology, Queen Square, London, United Kingdom

*For correspondence: yujiang.wang@newcastle.ac.uk

Abstract

Normative mapping is a framework used to map population-level features of health-related variables. It is widely used in neuroscience research, but the literature lacks established protocols in modalities that do not support healthy control measurements, such as intracranial EEG (icEEG). An icEEG normative map would allow researchers to learn about population-level brain activity and enable comparison of individual data against these norms to identify abnormalities. Currently, no standardised guide exists for transforming clinical data into a normative, regional icEEG map. Papers often cite different software and numerous articles to summarise the lengthy method, making it laborious for other researchers to understand or apply the process. Our protocol seeks to remedy this gap by providing a dataflow guide and key decision points that summarise existing methods. This protocol is used heavily in published works from our own lab (twelve peer-reviewed journal publications). Briefly, we take as input, icEEG recordings and neuroimaging data from people with epilepsy who are undergoing evaluation for resective surgery. As final outputs, we obtain a normative icEEG map, comprising signal properties localised to brain regions. Optionally, we can also process new subjects through the same pipeline and obtain their z-scores (or centiles) in each brain region, for abnormality detection and localisation. To date, a single, cohesive, dataflow pipeline for generating normative icEEG maps, along with abnormality mapping, has not been created. We envisage that this dataflow guide will not only increase understanding and application of normative mapping methods, but will also improve the consistency and quality of studies in the field.

Key features

- Resultant normative maps can be used to test a broad range of hypotheses in the neuroscience field
- Provides a more detailed walkthrough of the methods in the normative mapping study conducted by Taylor *et al.* (1) and other related publications from our lab (2–12)
- Offers flexibility – readers can tailor the final output by considering key Decision Points included throughout the protocol
- Involves sub-pipelines which may be useful to researchers in isolation (i.e., icEEG electrode localisation and/or interictal segment selection)

Keywords: Data pipeline, computational methods, neuroscience, normative mapping, abnormality mapping, neuroimaging, EEG, epilepsy


This protocol is used in: Brain (2022), DOI: 10.1093/brain/awab380; Journal of Neuroscience Methods (2024), DOI: 10.1016/j.jneumeth.2024.110180; eBioMedicine (2023), DOI: 10.1016/j.ebiom.2023.104848; Brain Communications (2023), DOI: 10.1093/braincomms/fcad292; Scientific Reports (2023), DOI: 10.1038/s41598-023-39700-7; Epilepsia (2023), DOI: 10.1111/epi.17663; Epilepsy & Behaviour Reports (2024), DOI: 10.1016/j.ebr.2024.100715.

The spatial localisation and database pipelines of this protocol are used in: Epilepsia (2020), DOI: 10.1111/epi.16580; NeuroImage: Clinical (2018), DOI: 10.1016/j.nicl.2018.01.028; Nature Communications (2024), DOI: 10.1038/s41467-024-52769-6; Epilepsia (2023), DOI: 10.1111/epi.17503; Epilepsia (2024), DOI: 10.1111/epi.18061.

Graphical overview

Starting materials per subject	
Required	Desirable
icEEG recordings	Post-surgery T1w MRI
Pre-implantation T1w MRI	Outcome data
Post-implantation CT scan	Channel details (spiking etc.)
Basic implantation details	More detailed subject metadata

Flowchart of pipelines involved in this protocol. Temporal Processing (Pipeline A) and Spatial Processing (Pipeline B) can be run sequentially or in parallel.

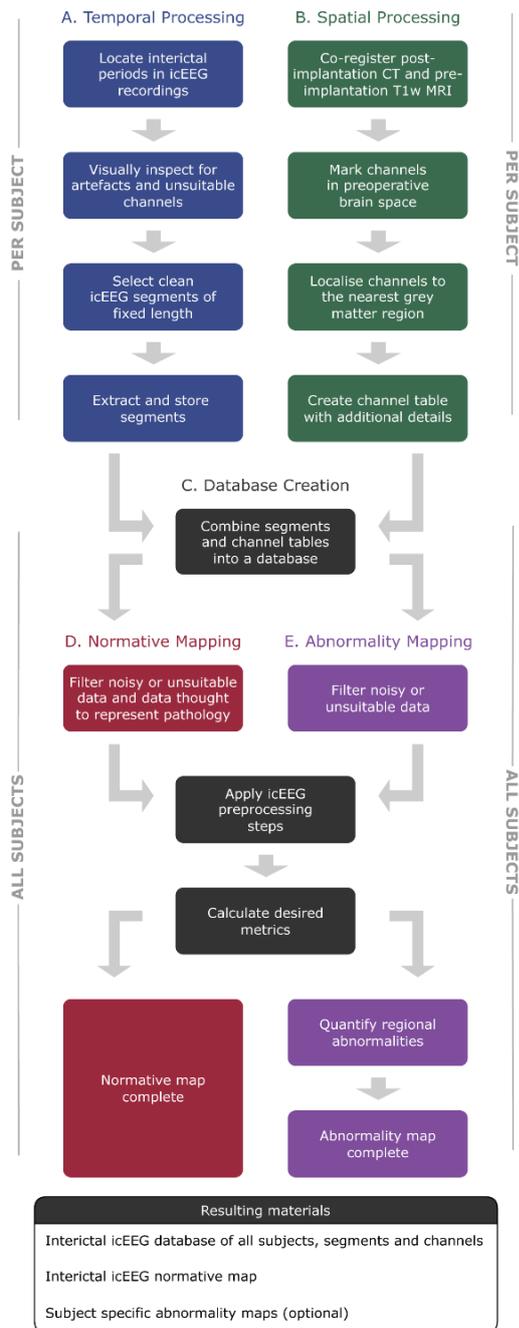

Acronyms

- **CT scan:** Computed tomography scan
- **EEG:** Electroencephalogram
- **GUI:** Graphical user interface
- **icEEG:** Intracranial EEG
- **MI:** Mutual information
- **PSD:** Power spectral density
- **RBP:** Relative band power
- **ROI:** Region of interest (usually in cortical grey matter or amygdala/hippocampus in this work)
- **SOZ:** Seizure onset zone
- **T1w MRI:** T1-weighted magnetic resonance imaging

Glossary

- **Channel:** One of several contact points found on an electrode, which records electrical brain activity.
- **Electrode:** A device comprising several recording channels, which is surgically implanted – types of electrodes include depth and grid.
- **Harmonics:** Power line noise which occurs at multiples of the initial frequency. For example, if there is power line noise at 50 Hz, it may also be present at 100 Hz and 150 Hz.
- **Period (of data):** A subsection of icEEG data taken from a subject's full dataset, which is taken forward for further investigation/possible segment selection.
- **Pre/post-implantation:** Refers to before/after the implantation of intracranial EEG electrodes. The word operation is avoided to prevent confusion between the electrode implantation operation, and the epilepsy surgery operation.
- **Pre/post-surgery:** Refers to before/after any resective epilepsy surgery. The word operation is avoided to prevent confusion between the electrode implantation operation, and the epilepsy surgery operation.
- **Resected (ROIs):** ROIs which were removed during epilepsy surgery.
- **Segment (of data):** A clean, interictal icEEG segment of fixed length that is inputted into the database for use in normative (or abnormality) mapping.
- **Spared (ROIs):** ROIs which were not removed during epilepsy surgery.
- **Subject:** An individual for which we have the necessary icEEG recordings and neuroimaging to proceed with the protocol. In our published works, this is primarily individuals with drug-resistant epilepsy undergoing evaluation for resective epilepsy surgery.

Background

Normative mapping is a framework in which features of health-related variables are mapped at the population level. A common example is growth charts, used to evaluate whether a child's growth falls within the normal range for their age. In neuroscience research, this framework is widely used to map normative ranges of brain activity, for example.

However, normative mapping is particularly challenging in modalities that do not support measurements in healthy control subjects, such as the highly invasive intracranial EEG (icEEG). An icEEG normative map of brain activity is valuable because it allows us to better understand populations of interest. Further, individuals can be compared to the normative map, helping to identify, characterise and localise potentially pathological abnormalities. Such work is available in the literature (1–7,13–20), revealing (for example) that icEEG normative mapping can localise epileptogenic tissue (1), and that the abnormality maps which follow are temporally stable (5).

Despite these promising findings, there is currently no detailed guide available for transforming clinical data, such as neuroimaging and icEEG recordings, into a normative map of brain activity localised to standardised brain regions. This literature gap hinders the consistent application of the methodology across studies. To remedy this, our protocol consolidates existing methods and provides key decision points for constructing normative maps comprising regional signal properties, along with (optionally) regional abnormality maps for new subjects.

One of the protocol's strengths lies in its modular design. Distinct sub-pipelines, such as the spatial processing pipeline, offer standalone utility. For instance, our own lab has employed this pipeline to demonstrate that the incomplete resection of the icEEG seizure onset zone (SOZ) is not associated with post-surgical outcomes (10), amongst other results (8,9,11,12). Furthermore, this protocol is written as a dataflow guide, allowing for flexibility. Researchers can choose their preferred programming language and data management software for implementation; decide whether to compute abnormalities; and apply the temporal processing pipeline to other modalities such as scalp EEG (for instance). Our protocol allows for adaptation to best suit the reader.

While the protocol provides a standardised framework, certain aspects, such as manual resection mask delineation, could benefit from modernisation. Automated methods exist within our own lab (21) and externally (22–25), although their reliability varies (26). The flexible, guide-like nature of our protocol is both a strength and a weakness. Importantly, we offer pseudo code and emphasize dataflow over low-level implementation, allowing researchers to tailor the process to their needs while fostering a deeper understanding of the method and related literature.

This protocol provides a cohesive dataflow on icEEG normative and abnormality mapping, which is currently missing from the literature. We aim to equip neuroscience researchers with the tools needed to develop normative maps, allowing them to explore new hypotheses consistently and effectively.

Equipment

1. PC or laptop capable of running the required software (below) with sufficient storage space. As a rule of thumb, we recommend at least 1GB per subject.

Software and datasets

Required

1. icEEG recordings for each subject
2. Post-implantation CT scan for each subject
3. Pre-implantation volumetric T1w MRI for each subject
4. Programming environment (e.g., MATLAB, Python, R)
5. Software to view icEEG data (e.g., EDFBrowser or the reader's chosen programming environment)
6. FreeSurfer (27) (we used version 7.3)
7. Electrode localisation tools such as `img_pip` (28)
8. Database software (e.g. MongoDB)

Optional/recommended

9. Post-surgery T1w MRI for subjects who proceeded to resective epilepsy surgery
10. ANTS toolbox (29)
11. 3D image viewer such as FSLEyes or ITK-Snap
12. RAMPS pipeline (21)

Note: If the reader is looking for suitable data, some of our published works (1–7) use the publicly available RAM dataset, found at <https://memory.psych.upenn.edu/RAM>.

Prerequisites

Readers intending to apply the protocol should have some experience working with neuroimaging and EEG data, along with some coding skills (in e.g., Python, MATLAB or similar).

Procedure

The folder structure throughout this protocol is the reader's preference. **Error! Reference source not found.** shows our recommendation for organising a subject's database and normative mapping files. Note, that it is common to have data from multiple batches, in this case often hospitals. Throughout the protocol, we will implement unique subject identifiers, a combination of a four-letter abbreviation of the hospital and the subject's index within that hospital. For example, the first subject from University College London Hospital would have the identifier `UCLH_001`. This will be referred to throughout as the Subject ID.

Our protocol comprises five sub-pipelines. Temporal Processing (Pipeline A) and Spatial Processing (Pipeline B) can be run sequentially or in parallel. Both pipelines are repeated for each subject. The remainder of the protocol (Pipelines C-E) is run with all subjects simultaneously. The outputs of Pipeline A & B are combined to create a database (Pipeline C). The database can be queried for subjects suitable for Normative Mapping (Pipeline D). Abnormality Mapping (Pipeline E) is optional and typically applied to either new subjects or a subset of existing subjects which are held out following the stage of Database Creation (Pipeline C).

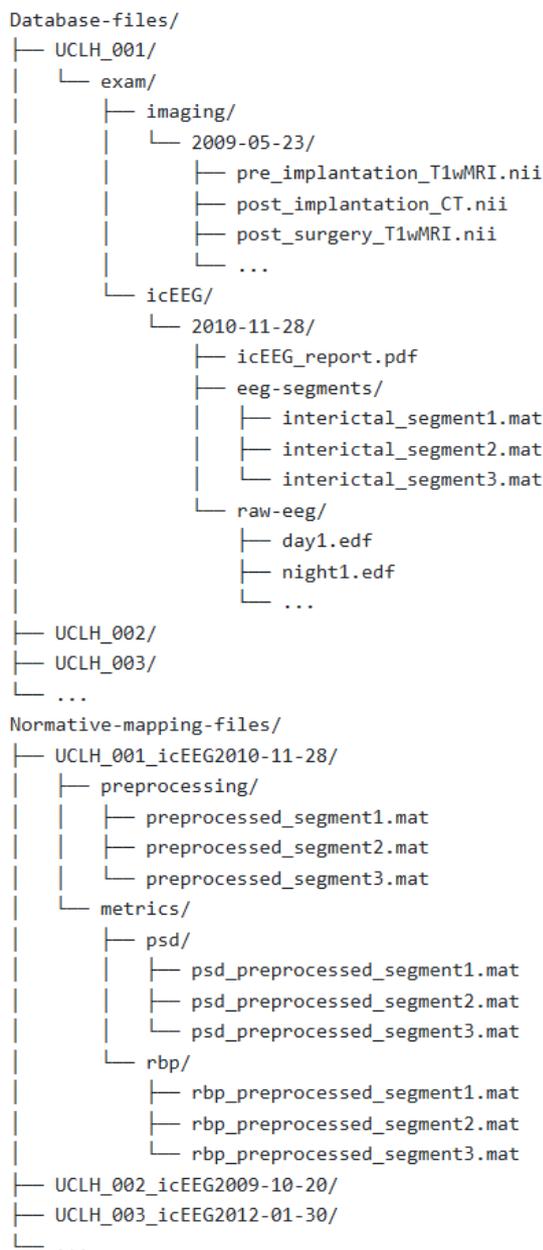


Figure 1: Example folder structure for Database Creation (Pipeline C) and subsequent Normative Mapping (Pipeline D)

A. Temporal Processing – icEEG segment selection

The **goal of Pipeline A** is to locate, extract, and systematically organise and store interictal segments from each subject. Before starting this pipeline, some **Decision Points (DPs)** must be addressed. These are the reader's choice and will affect the resultant normative map. For each, we provide the decisions made in Taylor *et al.* (1) as an example.

- **Decision Point A1:** What state of consciousness is the normative map? Taylor *et al.* (1) construct a normative map of relaxed wakefulness.
- **Decision Point A2:** What constraints will be imposed regarding seizures to ensure selected data periods and segments are truly interictal? In Taylor *et al.* (1), interictal data must be a minimum of two hours from any seizure activity.
- **Decision Point A3:** Which sampling frequency will the interictal segments use to create the normative map? Taylor *et al.* (1) ensure all segments are at least sampled at 200+ Hz. Data with a higher sampling frequency than the intended level can still be used and downsampled at a later stage. However, icEEG data with a lower sampling frequency than the intended level should be excluded.
- **Decision Point A4:** What is the desired referencing system? Taylor *et al.* (1) use common average. This can still be changed later if data is saved in referential montage (as recorded), but the reader should have a default option.
- **Decision Point A5:** What is the duration of the interictal segments used to construct the normative map? Segment length ideally would be consistent across all segments and subjects. Taylor *et al.* (1) uses 70-second segments. However, the reader could also choose to use segments of differing length if it is absolutely not avoidable, and co-vary for the length later statistically.
- **Decision Point A6:** How many interictal segments will be selected per subject? This will depend on how much data is available and may not be attainable for every subject. Multiple segments are selected to provide back-ups if a segment is excluded during preprocessing and can also be used for later validation. Taylor *et al.* (1) selected three segments per subject where possible.
- **Decision Point A7:** Will there be any time separation constraints on the selection of multiple segments? Taylor *et al.* (1) ensured that segments selected from the same subject and icEEG exam were separated by a minimum of four hours.

Once **Decision Points A1-7** have been finalised, begin Temporal Processing (Pipeline A). This is done on a *per subject* basis and assumes subjects have labelled, long-term icEEG recordings comprising various wake states and periods of ictal and interictal activity. If the starting data is more tailored for purpose, e.g., it is all interictal, some of steps 1-4 may

be skipped. If at any step, there is no data which meets the criteria, exclude this subject from further processing and begin with a new subject.

For a given subject...

1. Locate all periods within the icEEG files where the subject is in the required state of consciousness (**DP A1**) using icEEG file labelling and annotations.
Quality control: A time constraint can be imposed to be more rigorous with the state of consciousness criteria. For example, only locate periods between e.g., 8 am – 10 pm when looking for awake recordings. Alternatively, apply a data-driven method of detecting sleep stages (30,31).
2. Discard any periods where the interictal constraints are not met (**DP A2**) using icEEG file labelling and annotations. **Error! Reference source not found.** provides an example timeline of an icEEG exam for one subject. Following certain DPs being made, the figure demonstrates where suitable interictal periods would be located. **Error! Reference source not found.**A provides an example of an ictal period which must be avoided due to proximity to a seizure.
3. Discard any periods where the minimum sampling frequency is not met (**DP A3**).
4. Discard any periods which are shorter than the duration of the chosen segment length (**DP A5**).
5. Visually inspect the interictal periods for artefacts and noise, and discard any noisy periods. **Error! Reference source not found.**B demonstrates non-physiological noise which may be present in some icEEG files. This can arise, for example, when channels temporarily disconnect and should be avoided.
6. Retain only clean periods during the desired consciousness state. **Error! Reference source not found.**C provides an example of a suitable, clean interictal period which should be carried forward for further processing.

Quality control: If the subject has enough other periods available at this stage, avoid using the first 24 hours (for example) post-implantation. This will minimise immediate implantation and anaesthesia effects. See **Error!**

Reference source not found..

Troubleshooting: If data looks consistently noisy for a subject, before excluding this subject we suggest checking:

- Plotted amplitude is appropriate
- Necessary filters have been applied
- Desired re-referencing has been applied for visual inspection (**DP A4**)
- Power line noise has been removed (note, power line noise can differ by country)

If the data still looks noisy, it likely needs to be excluded.

Pause point: At this stage, the reader should have identified periods of clean, interictal data for the subject. How many periods are available will vary per subject, as it is impacted by factors such as the number of seizures they had during the icEEG examination, and how long the examination was.

7. Select the final segments of fixed length (**DP A5**, **DP A6**) from the periods of interictal data identified for this subject, ensuring all the imposed constraints are met (**DP A7**). Readers should check the power spectral densities (PSDs) as well as the icEEG when choosing final segments. **Error! Reference source not found.** demonstrates where segments might be selected from within suitable periods.

Quality control: If a segment has a small portion of ‘unsuitable’ channels, but is otherwise usable, make a note of the unsuitable channels and they can be removed further downstream. Unsuitable channels include those which are

noisy, faulty or showing non-biological recordings. Such channels can be identified through visual inspection of icEEG (see **Error! Reference source not found.D**), or through visual inspections of PSDs (see **Error! Reference source not found.E**). PSDs should follow a 1/f shape and may have a narrow spike indicating power line noise which can be dealt with at a later stage. Channels exhibiting a different pattern should be marked as unsuitable. More rigorous checks are carried out later, but visual checks are still useful here. **Error! Reference source not found.F** shows a clean, interictal segment.

Quality control: Some channels in the interictal segments may exhibit epileptic activity, e.g. spikes. If any spiking channels are identified during segment selection, either visually or using a spike detection algorithm (32,33), they should be noted as such so they can be excluded further downstream.

Pause point: At this stage, the reader should have identified a selection of clean, interictal segments for the subject.

8. Create a suitable folder to store this subject's extracted segments. See **Error! Reference source not found.** for our recommendation.
9. For one segment use your preferred programming language to extract the following data from the original file format (e.g., extract the relevant data from a .edf file and store it as a .mat file). These should be stored within the folder created in the previous step.
 - a) All channel names
 - b) Names of channels labelled as unsuitable (a subset of a))
 - c) Names of channels labelled as spiking (a subset of a))
 - d) Time series data (a matrix of size channels \times time points)
 - e) Date and time at chosen precision (e.g. second) at the start of the segment
 - f) Duration of the segment
 - g) Sampling frequency
 - h) State of consciousness
 - i) Subject ID
 - j) Hospital at which the data was recorded
 - k) Segment number
10. Repeat step 9 for all remaining segments.

Repeat the 10 steps above for all subjects.

Endpoint: A selection of interictal, fixed length icEEG segments in the desired state of consciousness extracted for each subject, with metadata included and unsuitable channels noted.

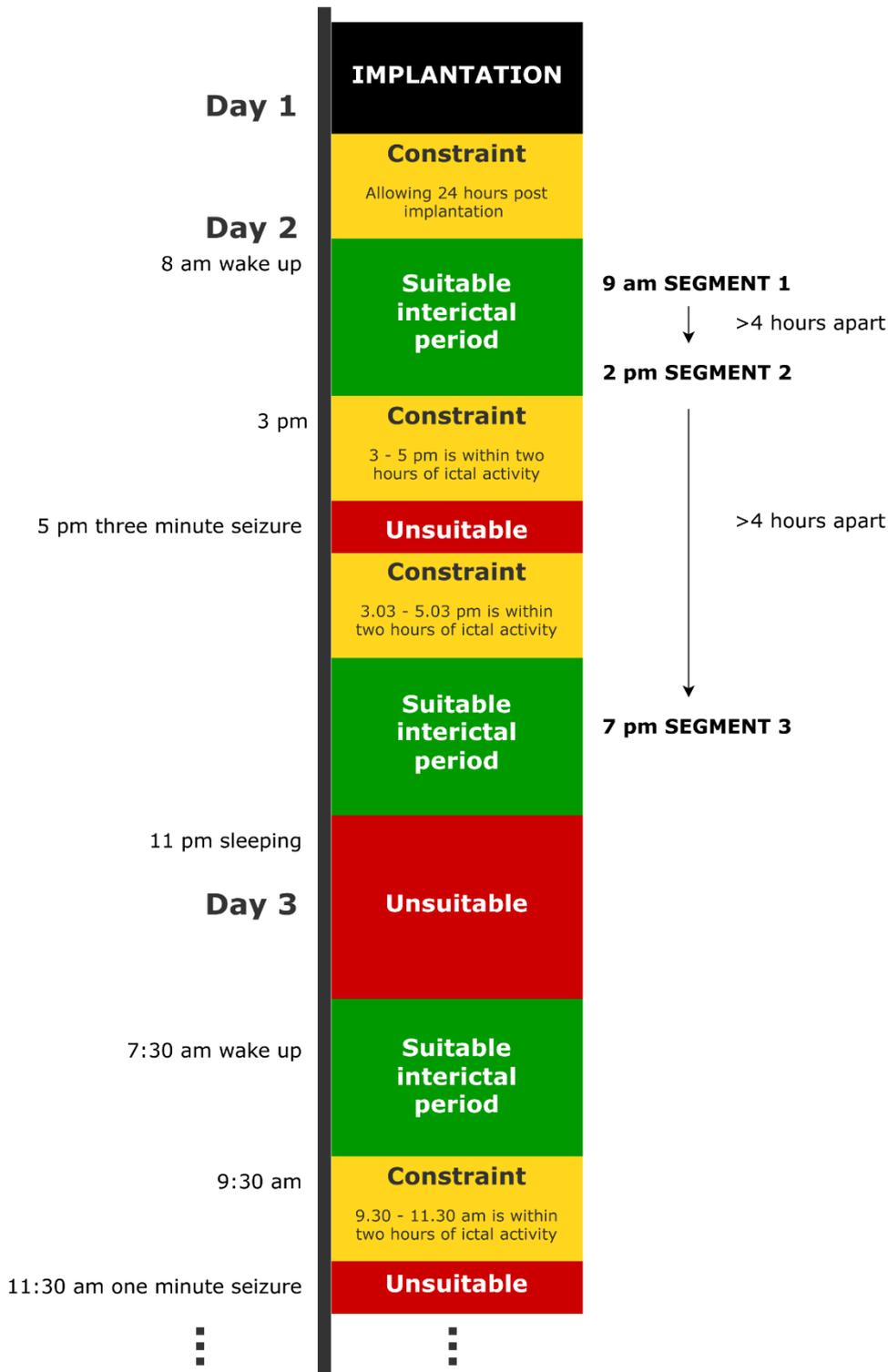


Figure 2: An example timeline of an icEEG exam for a dummy subject, demonstrating where suitable periods of interictal data would be, and the temporal location of the final interictal segments. In this example, the normative map is being constructed
 Bio-protocol Manuscript Template, Version 5 (June 27, 2023)

in the wake state (DP A1) and three segments are being selected (DP A6) under the constraint that they are 2 hours away from seizures (DP A2), 4 hours away from one another (DP A7) and over 24 hours after implantation (Quality control).

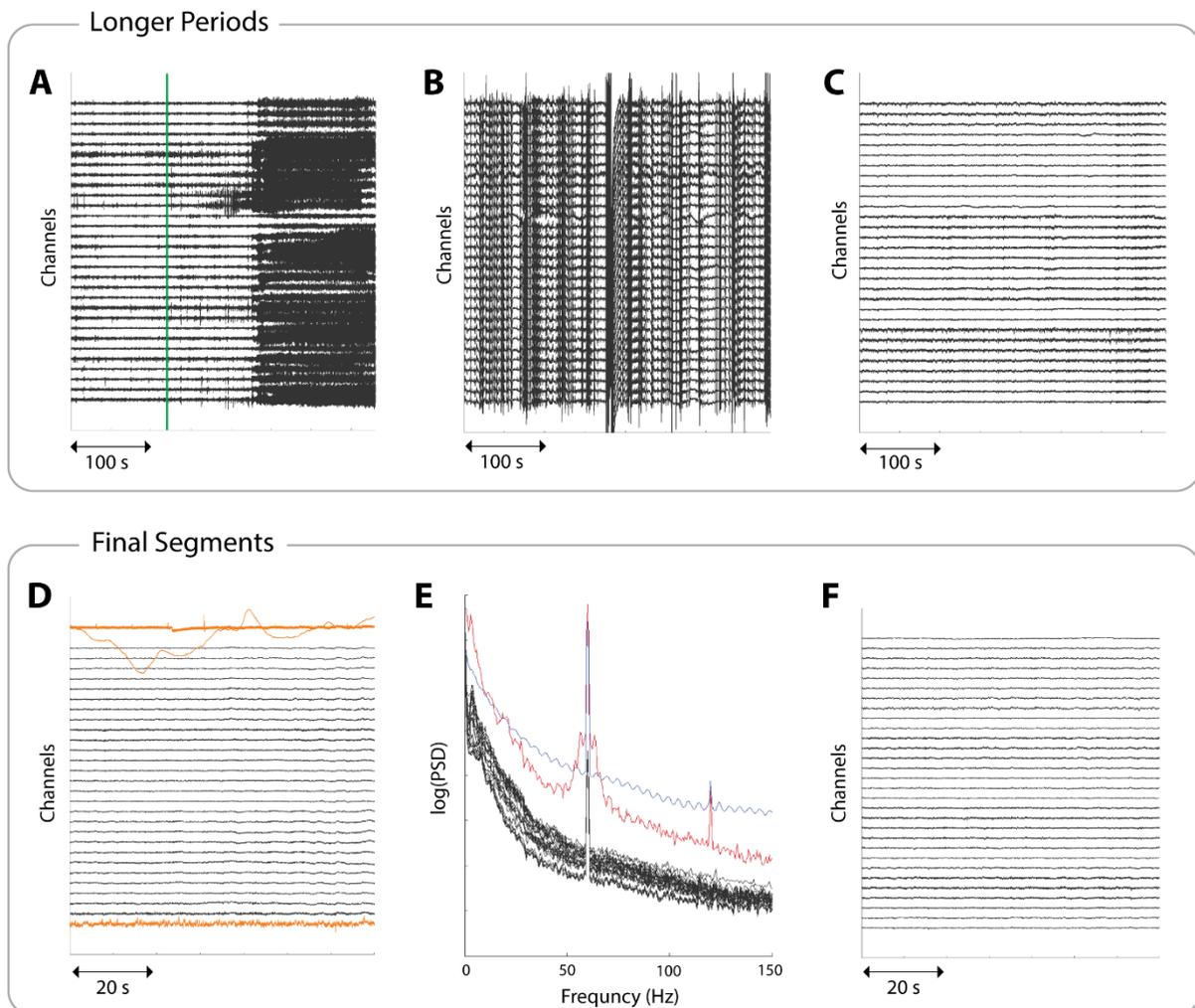


Figure 3: Example icEEG data to aid in identifying suitable periods and segments. Each panel shows a different subject and a subset of 30 channels. **A)** A 380 s data period involving a seizure. The seizure onset, as marked by clinicians, is indicated by a green line. This line would correspond to the start of a red 'Unsuitable' box in Figure 2. **B)** A 380 s data period showing non-physiological noise, indicated by the 'streaky' appearance. Such large amplitudes are not representative of brain activity and data exhibiting this should be discarded. **C)** A clean 380 s data period, ready to proceed to segment selection. This would correspond to 'Suitable interictal period' boxes in Figure 2. **D)** A 70 s interictal segment with three examples of unsuitable channels highlighted in yellow. The segment is usable as long as these channels are noted so they can be removed downstream. **E)** The power spectral densities (PSDs) of a 70 s interictal segment. Suitable channels are coloured black, showing a $1/f$ curve with power line noise at 60 Hz. The red line shows a channel with a broad peak around 60Hz and an overall higher power in most frequencies (unlike the black channels), which most likely indicates a faulty channel. The blue line shows a channel which likely represents non-physiological activity due to its oscillatory pattern in the frequency domain. Both channels should be marked as unsuitable. **F)** A clean 70 s interictal segment ready to be extracted and stored.

B. Spatial processing – electrode localisation and resection mask delineation

The **goal of Pipeline B** is to localise electrodes to standard brain regions of interest (ROIs), record additional channel details and, where applicable, create resection masks. Before starting this pipeline, some **Decision Points (DPs)** must be addressed. These are the reader's choice and will affect the resulting normative map. For each, we provide the decisions made in Taylor *et al.* (1) as an example.

- **Decision Point B1:** What constraint will be imposed for excluding channels based on distance from grey matter? Taylor *et al.* (1) exclude channels which appear to be located more than 5mm from grey matter.
- **Decision Point B2 (optional, requires resection mask):** What threshold will be implemented to define a channel as recording from subsequently resected tissue? Taylor *et al.* (1) define channels as recording from resected tissue if they were within 5mm of the resection mask.

Once **Decision Points B1 & B2** have been finalised, begin Spatial Processing (Pipeline B). This is done on a *per-subject* basis and assumes subjects have a pre-implantation T1w MRI, a post-implantation CT scan and, if drawing resection masks, a post-surgery T1w MRI. Clinical reports may have information on which channels were recording from subsequently resected tissue. If so, **Decision Point B2**, step 5 and steps 7.5-7 can be considered unnecessary. Be aware that not all subjects proceed to resective epilepsy surgery. Subjects without resection masks or resection details may still be useful, depending on the subsequent research questions. If at any step, there is no data which meets the criteria, exclude this subject from further processing and begin with a new subject.

For a given subject...

1. Run the FreeSurfer 'recon-all' pipeline on the pre-implantation 3D T1-weighted MRI.

Note: Following successful completion, the file `mri/orig.mgz` will be generated. This file has 1mm isotropic voxels and is the reference space for all future analyses. This step also generates the file `mri/aparc+aseg.mgz`, which contains the volumetric parcellation of different ROIs. This uses the Desikan-Killiany parcellation (34) built into FreeSurfer, and is used later in step 7.

Quality control: Check the outputs of step 1 and consider performing manual edits via control points where appropriate (35).

2. Co-register the post-implantation CT scan to the `orig.mgz` file. This may require conversion from `mgz` format to `nifti (nii.gz)` format. For this, the FreeSurfer `mri_convert` command can be used. Registration should be rigid body (6 degrees of freedom), without deformation. The ANTS toolbox (29) contains a highly effective set of tools for this purpose. The following command can be used from the toolbox:

```
/usr/local/ANTs/bin/antsRegistration --verbose 1 --dimensionality 3 -
-float 0 --collapse-output-transforms 1 --output [
<Output_FilePaths>] --interpolation Linear --use-histogram-matching 0
--winsorize-image-intensities [ 0.005,0.995 ] --initial-moving-
transform [ <input_orig.nii.gz>,<input_CT.nii.gz>, 1 ] --transform
Rigid[ 0.1 ] --metric MI[<input_orig.nii.gz>,
<input_CT.nii.gz>,1,32,Regular,0.25 ] --convergence [
1000x500x250x100,1e-6,10 ] --shrink-factors 8x4x2x1 --smoothing-
sigmas 3x2x1x0vox
```

The rigid registration uses mutual information (MI) to measure concordance between the images. Since the images are from different modalities, contrast for specific tissues will be altered. MI compares histograms of the two images and aligns them, detecting similar anatomical patterns between the images, even if they have contrasting intensities. This method is widely used for multi-modal registration.

Quality control: Following registration, visually inspect that the images are correctly aligned using software for viewing 3D imaging data such as FSLeyes or ITK-Snap (see Figure 4A for an example of a correctly aligned image).

3. Manually mark channels in pre-operative brain space using, for example, the open-source Python package `img_pipe` (28).

Note: A full protocol for this software has been published previously (28). Here, we summarise the section “Electrode Identification on the Co-registered CT Scan” to obtain Cartesian coordinates for each channel.

1. Follow the instructions in Hamilton *et al.* (28) to launch an interactive Python GUI showing the registered CT scan overlaid on the skull-stripped MRI.
2. Navigate the crosshairs to channel 1 of an electrode in the brain. Channel 1 is usually the deepest point of contact, but double-check implantation information to confirm this.
3. Add a new electrode by pressing the ‘n’ key, triggering a prompt to name it. For simplicity, we recommend following the clinician’s naming systems. This is often anatomical abbreviations (HIPPI, AMYG...) or simple alphabetical labelling (A, B...).
4. Press the ‘e’ key to add a channel at the crosshair position. A coloured shape will appear in the GUI indicating where the channel has been marked, along with a legend showing the electrode name. The channel coordinates will be automatically saved in a folder named after the electrode.

Quality control: Ensure the crosshairs are in the centre of the channel in *all* views in the GUI (coronal, sagittal, axial).

5. Continue to mark the channels on the current electrode by pressing ‘e’ over each one. Work in order up the electrode from the deepest to the most superficial channel.
 6. Repeat steps 3.2-3.5 above for each electrode in the brain.
4. Once the channels are marked for all electrodes in the subject’s brain, export the channel locations. If using the pipeline in Hamilton *et al.* (28), these will be in the same space as the `orig.mgz` file described in step 1.

Quality control: Write custom code or follow the guidelines in Hamilton *et al.* (28) to produce a rendering of the electrodes in the subject’s brain space. Visually confirm that the electrodes are in the expected locations by comparison to implantation information, Figure 4B provides an example of this.

Pause point: At this stage, the reader should have a folder containing channel locations for each of the subject’s implanted electrodes.

5. [optional] Generate a resection mask in pre-implantation `orig.mgz` space. We recommend the RAMPS pipeline (21). In brief, this resection mask pipeline generates a mask by ...
 1. Registering a post-surgery T1w MRI to the pre-implantation T1w MRI (using the ANTS toolbox (29))
 2. Performing tissue segmentation to identify a resection cavity (using ANTS ATROPOS methods (29))
 3. Performing multiple rounds of segmentation, dilation and erosion to generate an accurate delineation

Once installed, the RAMPS pipeline can be executed with the following command:

```
python /Path_to/RAMP.py </Path_to/Pre-implantation-Scan.nii.gz>
</Path_to/Post-surgery-Scan.nii.gz>
</Path_to_Output_Folder_file_path/> <Output_Prefix> <Hemisphere>
<Lobe>
```

where `<Hemisphere>` is the hemisphere containing the resection (e.g. “L” or “R”), and `<Lobe>` is the lobe(s) of resection (e.g. “T” for temporal, or “F” for frontal).

Quality control: Overlay the resection mask with the `orig.mgz` file to ensure it is in the expected location. Figure 4A includes a visualisation of a resection mask.

6. Generate a table for each channel which contains the following information:
 - a) Name of the channel (electrode name, followed by channel number, e.g. HIPPI1)
 - b) Electrode type (such as “G” for grid electrodes, “D” for depth electrodes)
 - c) Hemisphere from which each channel records from: (“Left”, “Right”)

7. Write code in the reader's preferred programming language (e.g., MATLAB) to:
 1. Load the channel locations (generated in step 3) and table containing channel information (generated in step 6).
 2. Load the `aparc+aseg.mgz` file (generated in step 1), which contains volumetric grey matter ROIs, and a volumetric segmentation of the white matter. A surface representation of the parcellation is shown in Figure 4C.

Note: In MATLAB the FreeSurfer function `MRIread.m` can be used for this.
 3. Localise channels to the nearest grey matter region in the same hemisphere by minimising the Euclidean distance.

Note: It is possible that a channel implanted between the two hemispheres may be closer to a region in the opposite hemisphere than that from which it records. It is therefore important to constrain by known side.
 4. Exclude channels that are more than the specified distance from grey matter (**DP B1**)/are located in white matter by setting their assigned region to 'n/a'.
 5. [optional] Load the resection mask generated in step 5.
 6. [optional continued] Compute the Euclidean distance from each channel to each voxel located in the mask.
 7. [optional continued] Save the minimum distance computed in the previous step.

Pause point: At this point, the reader should have localised channels to ROIs, excluding those in white matter or located too far away from grey matter. Optionally, for each channel, the reader should have recorded the minimum distance of each channel to the resection mask.

8. Load the output from steps 6 (channel coordinates & names) and 7 (ROI localisations and distances from channels to resection masks) into your preferred programming software.
9. Create a final table encompassing all channel details for the subject, with the variables/columns outlined below. Table 1 provides an example channel table for one electrode. Table 2 shows the desirable (but not essential) variables for the same example electrode.
 - a) Channel name (e.g., HIPP1, first recorded in step 3)
 - b) Coordinates (x, y, z) of the channel (first recorded in step 3) in the same space as `orig.mgz`
 - c) Type of electrode that the channel is on (step 6)
 - d) Hemisphere from which the channel records from (step 6, "Left", "Right")
 - e) Numerical index of the ROI to which the channel has been assigned using the Desikan-Killiany atlas. A mapping from numerical index to region is available in Section 2 of an ROI tutorial on the FreeSurfer Wiki (36). A visualisation of the atlas is provided in Figure 4C.
 - f) [optional] Binary indicator of whether the channel was recording from tissue that was subsequently resected (using the output from step 7.7) (**DP B2**)

The following variables are desirable, but not essential:

 - g) Binary indicator of whether the channel exhibited spikes in the interictal icEEG recordings (found in clinical reports if available, or may have been identified in Pipeline A)
 - h) Binary indicator of whether the channel is within any structural abnormality, such as a lesion (found in clinical reports if available, or identified through visual inspection of neuroimaging data)
 - i) Binary indicator of whether the channel is located within the clinician-defined SOZ

10. Save the resulting channel table for this subject using the subject identifier in the file name.

Note: With the exception of the original neuroimaging files and the final channel table, the reader is free to discard all other tables, variables and files created during this pipeline. The information in the final channel tables will ultimately be recorded in the database (created in Pipeline C). Hence, channel tables are not included in the example folder structure in **Error! Reference source not found.**

Repeat the 10 steps above for all subjects.

Endpoint: A channel details table for every subject, containing the name, location and regional localisation of every channel the subject had implanted, along with other channel-level metadata where available.

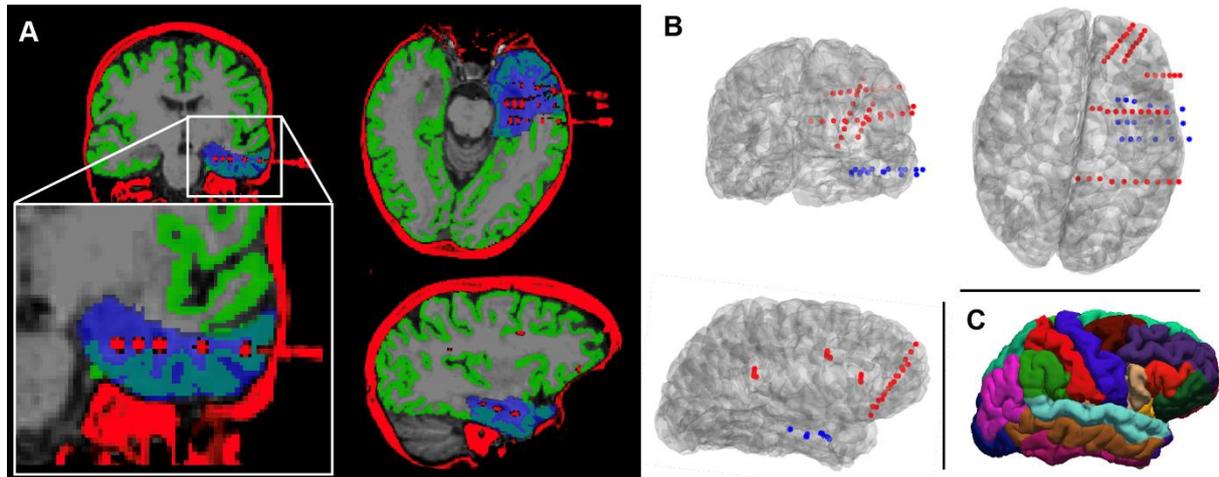


Figure 4: *Quality control checks and outputs from Pipeline B for an example subject. A) The correctly aligned post-implantation CT scan and orig.mgz file following co-registration (Quality control, step 2). Overlaid in blue is the resection mask, which delineates the tissue that was subsequently resected during surgery (Quality control, step 5). The mask was generated using the RAMPS pipeline (21). B) A rendering of the subject's implanted electrodes following channel localisation (Quality control, step 4). The channels on each electrode are indicated by coloured dots; blue channels were defined as recording from subsequently resected tissue, whilst red channels were spared (step 9). Looking at both panels A) and B), it can be determined that the subject had surgery on the temporal lobe. C) A visualisation of the Desikan-Killiany parcellation (34) to which channels will be localised.*

name	location orig			electrode type	hemisphere	ROI ID	ROI name
LHIPP1	-29.83	-15.10	-9.69	depth	left	17	Left-Hippocampus
LHIPP2	-34.91	-15.10	-8.60	depth	left	17	Left-Hippocampus
LHIPP3	-39.70	-15.10	-7.65	depth	left	17	Left-Hippocampus
LHIPP4	-44.59	-14.44	-6.77	depth	left	39	lh-middle-temporal
LHIPP5	-49.29	-14.44	-5.68	depth	left	39	lh-middle-temporal
LHIPP6	-54.36	-13.71	-4.86	depth	left	39	lh-middle-temporal
LHIPP7	-58.97	-13.71	-4.29	depth	left	39	lh-middle-temporal
LHIPP8	-64.14	-13.24	-3.15	depth	left	39	lh-middle-temporal
LHIPP9	-68.00	-13.24	-2.46	depth	left	39	lh-middle-temporal
LHIPP10	-73.33	-13.24	-1.29	depth	left	39	lh-middle-temporal

Table 1: A channel table for one electrode targeting the hippocampus in an example subject.

name	Channel recorded from tissue that was...			
	subsequently resected	exhibiting spikes	structurally abnormal	within the clinician-defined SOZ
LHIPP1	0	1	0	1
LHIPP2	0	1	0	1
LHIPP3	0	1	0	0
LHIPP4	0	1	0	0
LHIPP5	0	1	0	0
LHIPP6	0	0	0	0
LHIPP7	0	0	0	0
LHIPP8	0	0	0	0
LHIPP9	0	0	0	0
LHIPP10	0	0	0	0

Table 2: A continuation of Table 1, showing an example of the additional, desirable variables that should be recorded in the channel tables if available.

C. Database creation

The **goal of Pipeline C** is to create an organised database which can be queried for downstream work.

Database creation aggregates *all subjects' data* and uses the outputs from Temporal Processing and Spatial Processing (Pipelines A & B). This database will include all subjects' data and channel information (including channels marked as unsuitable, resected channels, etc). The choice of database software is the reader's personal preference. In our case, we have used MongoDB as a NoSQL database for increased flexibility and scalability in a range of research applications, not limited to normative mapping. In theory, however, for the purposes of this paper, relational databases, or even using several spreadsheets could also be a suitable method of databasing. The user may wish to use relational databases if the database is limited in its purpose (e.g. only for normative mapping), and a stronger guarantee of data normalisation is required. We would not recommend using several spreadsheets if the database is highly complex with one-to-many or many-to-many relationships, e.g. where several subjects have repeat visits or second surgeries, or when several hospital sites and subjects are involved.

Instead of a step-by-step process, we here describe the structure of our database in detail, to allow readers to create similar structures using their preferred software. See Figure 5 for a visualisation of the database structure. At the end of this section, we also include the Object-Data-Map for MongoDB (written with Python and mongoengine syntax) for completeness and clarity.

To begin, we list some basic database requirements, which apply regardless of software choice:

- The database consists of various *classes* which define its structure
- Classes have properties, which we call *fields*, where we input information
- Fields can be *required*, meaning we must have input to make an entry, or they can be *optional*
- To ensure data quality, field inputs can be subject to certain *constraints*, meaning they can only take certain values

Specifically, we have outlined the structure for our icEEG database below. As in Figure 5, required fields are **bold**, and constraints are indicated by *italics*. Additionally, the **font is changed** to indicate database variables within the text. Coloured text indicates the class name and is level specific.

Level One

Level One: Subject

Subject is the highest-level class. We must have the **Subject** class to be able to input any lower-level data, such as the icEEG segment information. All lower-level classes are optional, we can simply record knowledge of an individual subject without any further specifics on their examinations or treatments.

Required fields in the **Subject** class are **Hospital** and **Subject ID**. Here **Hospital** denotes the site at which the subject was recorded; we choose to constrain it to a 4-letter key, for example, “ULCH” for “University College London Hospital”. As mentioned previously, **Subject ID** combines the individual’s **Hospital** code with a unique numerical identifier. The **Subject ID** appears in all subsequent classes, as such we can follow each subject throughout the levels of the database.

Optional fields are reader preference, we include **Date of birth** (or simply **Year at birth**), **Sex at birth**, **Handedness** and **Age of epilepsy onset**. **Sex at birth** is an example of a field which should be constrained, taking exactly one entry out of *male* and *female*. Similarly, **Handedness** should be constrained to *left*, *right*, and *ambidextrous*.

Level Two

Level Two comprises three classes which branch off from **Subject** (Level One). These are **Epilepsy Status**, **Exam** and **Treatment**. There is an optional, one-to-many relationship between **Subject** and each of these lower-level classes, as indicated by the connecting lines used in Figure 5. This means a given subject could have any number of entries for each Level Two class (including no entry). Each Level Two class is outlined in more detail below.

Level Two: Epilepsy Status

Epilepsy Status records information about the state of the individual's medical condition. This would for example reflect the epilepsy status at a particular point of follow-up with the subject, e.g. capturing if they had seizures in the recent past. A **Subject** may have multiple records of their epilepsy status across a few years, or they may have no record of epilepsy status at all.

Required fields in the **Epilepsy Status** class are **Subject ID** and **Date of status**. The **Date of status** field must be unique within **Subject ID**.

Optional fields we include are **Status epilepticus**, **Complex partial seizures** and **Secondary generalised seizures**, which all take a Boolean input of *true* if the individual had the seizure type, and *false* otherwise. Additional optional fields include the monthly frequency of each seizure type.

Level Two: Exam

Exam refers to the assessments performed, e.g. to try and localise the seizure origin. In our case, individuals often have both icEEG and neuroimaging exams available. Further, a **Subject** may have more than one icEEG recording session. Each different exam requires a new entry.

Required fields in the **Exam** class are **Subject ID**, **Exam ID**, **Exam type** and **Exam date**. Similarly to **Subject ID**, **Exam ID** is a string which is unique to each of the subject’s exams. It appears in all subsequent classes (Level Three) allowing the reader to trace back to the specific exam. **Exam type** is constrained to a list of values, such as *icEEG*, *MRI*

or *CT*.

For each **Exam type** there are further optional fields indicating the file location of the raw icEEG or neuroimaging files, along with any available related files such as the icEEG report. Note, that raw data is not stored in the database itself, just the location is recorded. For example, taking **Error! Reference source not found.** as a suitable folder structure, for subject UCLH_001's icEEG exam, the database entry for Raw icEEG folder location would be 'Database-files/UCLH_001/exam/icEEG/2010-11-28/raw-eeg'.

Level Two: Treatment

Treatment refers to the approach taken to manage epilepsy, such as surgery. Note, that treatment data might not be available for some individuals. Alternatively, an individual may have two resective epilepsy surgeries and thus will have two **Treatment** entries.

Required fields in the **Treatment** class are **Subject ID**, **Treatment type** and **Treatment date**. In our data, **Treatment type** is constrained to *Epilepsy surgery*, but in theory medication treatments can also be captured here. **Treatment date** must be unique for any combination of **Subject ID** and **Treatment type**.

Optional fields include **Outcome** fields (for ILAE, Engel or both) and **Outcome date**. In the case where **Treatment type** = *Epilepsy surgery*, there are additional optional fields including **Surgery type**, **Surgery pathology** and **Surgery hemisphere**. The latter is an example of a constrained field, taking one entry out of *left* and *right*.

Level Three

Level Three comprises three classes that branch off from **Exam** (Level Two) and only become available when the subject has an input for **Exam**, and that input includes **Exam type** = *icEEG*. The Level Three classes are **Channel**, **icEEG Segment** and **Seizure**. If an individual had more than one icEEG exam, the subsequent Level Three classes corresponding to each session will be distinguished using the **Exam ID** defined in Level Two.

There is an optional, one-to-many relationship between **Exam** and each of these lower-level classes, as indicated by the connecting lines used in Figure 5. This means that a particular subject's icEEG exam could have any number of entries for each Level Three class (including no entry). Each Level Three class is outlined in more detail below.

Level Three: Channel

Channel refers to the recording channels on each electrode the subject had implanted in the icEEG exam. Subjects can have any number of channels implanted, and each one has its own **Channel** input. This class is where all the information in the channel tables created in Spatial Processing (Pipeline B) is stored. For a given subject, each row of their channel table would correspond to one **Channel** entry.

Required fields in **Channel** are **Subject ID**, **Exam ID** and **Channel name**. The latter must be unique within this exam (using **Exam ID**) and is typically a combination of the electrode name and the channel's location on the electrode, e.g. in our data, HIPPI would be the deepest channel on the electrode targeting the hippocampus.

Optional fields include Channel hemisphere, Electrode type, Channel location, ROI ID, and ROI name. There are a few constraints on these fields:

- Channel hemisphere should be constrained to *left* and *right*
- Electrode type refers to the type of electrode the channel sits on and is constrained to *grid*, *strip* or *depth*
- The ROI is the brain parcellation that the channel is assigned to, and is identified using both a number (ROI ID) and a label (ROI name)
- ROI ID and ROI name should be one of those found in the Desikan-Killiany parcellation (34)

There is a further set of optional, Boolean fields referring to the suitability of the channel for normative mapping. These are indicators of whether the channel is Resected, Spiking, Within structurally abnormal tissue or Within the SOZ. Such information may be found in clinical reports or may be ascertained in previous pipelines.

Level Three: icEEG Segment

icEEG Segment stores information on the interictal segments selected in Pipeline A. Some subjects may not have any suitable interictal segments in their icEEG exam, whilst others may have had several selected (**DP A6**).

Required fields in **icEEG Segment** are **Subject ID**, **Exam ID**, and the **Date-time**, **Duration**, **Sampling frequency**, **Reference type**, **Wake status** and **File location** of the segment. In our data, the **Reference type** is constrained to *referential* and *bipolar*. The segment's **Wake status** should be constrained to the standard sleep stages (*Wake*, *N1*, *N2*, *N3*, *REM*), with *NA* as a final option for any segments where the state of consciousness could not be determined. Again, taking **Error! Reference source not found.** as a suitable folder structure, for subject UCLH_001's first interictal segment, the database entry for **Segment file location** would be 'Database-files/UCLH_001/exam/icEEG/2010-11-28/eeg-segments/interictal_segment1.mat'.

We do not include any optional fields in **icEEG Segment**.

Quality control: As an extra precaution, Decision Points from Temporal Processing (Pipeline A) can be quality control checked here. For example, are the inputted segments the desired duration (**DP A5**)? Are multiple segments separated by the outlined time constraints (**DP A7**)?

Level Three: Seizure

Seizure stores information on ictal events which may have occurred during the examination. Subjects may have any number of seizures during the icEEG exam (including none). Seizure data itself is not of interest for this pipeline, however knowing what type of seizures occur and when is crucial for identifying suitable interictal segments – any information that may be relevant to segment selection is therefore stored within this **Seizure** class.

Required fields in **Seizure** are **Subject ID**, **Exam ID** and **Seizure date-time**. The latter must be unique within the exam

(using **Exam ID**).

Optional **Seizure** fields include the seizure Duration, Type and Onset channels. Duration determines the endpoint of the seizure, helping to know when an interictal segment can start. Further, readers may want to impose different constraints on different seizure types. For example, perhaps interictal segments should be at least 2 hours away from subclinical seizures but must be at least 4 hours away from focal seizures (**DP A2**). Seizure type can help with this if the information is available. Finally, Seizure onset channels might be found in clinical reports, if so, this field can help identify channels which are unsuitable for normative maps.

That concludes the written description of the database, which can be queried for various subsets of data. Below, we include the Object-Data-Map we used for MongoDB as an example.

End point: An interpretable, organised database of all subjects which can be queried to identify segments based on various criteria, in particular, those suitable for normative mapping.

PIPELINE C: DATABASE CREATION

EXAMPLE DATABASE OBJECT DATA MAPPING

```
from mongoengine import *
```

INPUT OPTIONS FOR THE CONSTRAINED FIELDS

General Constraints

```
boolean = (('0', 'FALSE'),  
          ('1', 'TRUE'))
```

```
Hemispheres = (('L', 'Left hemisphere'),  
              ('R', 'Right hemisphere'))
```

Subject Constraints

```
Hospital_ID = (('UCLH', 'NHS UCL Hospitals'),  
              ('GOSH', 'Great Ormond Street'))
```

```
Sexes = (('M', 'born as male'),  
        ('F', 'born as female'))
```

```
Handednesses = (('L', 'left-handed'),  
                ('R', 'right-handed'),  
                ('A', 'ambidextrous'))
```

Exam Constraints

```
Exam_types = (('icEEG', 'Intracranial EEG examination'),  
             ('MRI', 'Magnetic resonance imaging'),  
             ('CT', 'Computer Tomography'))
```

Treatment Constraints

```
Treatment_types = (('Surgery', 'Epilepsy Surgery'))
```

Channel Constraints

```
electrode_types = (('grid', 'grid ECoG'),  
                  ('strip', 'strip ECoG'),  
                  ('depth', 'depth or SEEG'))
```

icEEG Segment Constraints

```
icEEG_reference_types = (('referential', 'as recorded relative to some reference'),  
                        ('bipolar', 'some bipolar montage'))
```

```
sleep_stages = (('W', 'awake'),  
               ('N1', 'non-REM light sleep, stage 1'),  
               ('N2', 'non-REM light sleep, stage 2'),  
               ('N3', 'non-REM deep sleep'),  
               ('R', 'REM sleep, vivid dreaming'))
```

Seizure Constraints

```
seizure_types = (('focal', 'Focal seizure'),  
                ('subclin', 'Subclinical seizure'),
```

```
('sg', 'Secondary generalised seizure'))
```

```
## SETTING UP CLASSES
```

```
## Level One: Subject Class
```

```
class subject(Document):
```

```
    # required fields
```

```
    Hospital = StringField(choices=Hospital_ID, required=(True))
```

```
    subjID = StringField(required=(True), unique_with=['Hospital'])
```

```
    # optional fields
```

```
    DOB = DateField()
```

```
    Sex_at_birth = StringField(choices=Sexes, max_length=1)
```

```
    Handedness = StringField(choices=Handednesses)
```

```
    Age_of_epilepsy_onset = FloatField()
```

```
    # link to Level Two classes
```

```
    treatments = ListField(ReferenceField('treatment'))
```

```
    exams = ListField(ReferenceField('exam'))
```

```
    epilepsy_status = ListField(ReferenceField('epilepsy_status'))
```

```
## Level Two: Epilepsy Status
```

```
class epilepsy_status(Document):
```

```
    # required fields
```

```
    subject_id = ReferenceField('subject', required=(True))
```

```
    date_of_status = DateField(required=(True), unique_with=['subject_id'])
```

```
    # optional fields
```

```
    status_epilepticus = BooleanField()
```

```
complex_partial_seizures = BooleanField()  
CPS_monthly_frequency = FloatField()
```

```
secondary_generalised_seizures = BooleanField()  
SGS_monthly_frequency = FloatField()
```

Level Two: Exams

```
class exam(Document):
```

```
    # required fields
```

```
    subject_id = ReferenceField('subject', required=(True))  
    exam_type = StringField(choices=Exam_types, required=(True))  
    exam_date = DateField(required=(True))
```

```
    # optional fields specific to exam_type = icEEG
```

```
    iceeg_report_file_location = StringField()  
    raw_iceeg_folder_location = StringField()
```

```
    # optional fields specific to exam_type = MRI
```

```
    pre_implantation_mri_file_location = StringField()  
    post_surgery_mri_file_location = StringField() # if available, for resection masks
```

```
    # optional fields specific to exam_type = CT
```

```
    post_implantation_ct_file_location = StringField()
```

```
    # link to Level Three classes
```

```
    channels = ListField(ReferenceField('channel'))  
    iceeg_segment = ListField(ReferenceField('iceeg_segment'))  
    seizures = ListField(ReferenceField('seizure'))
```

Level Two: Treatments

```
class treatment(Document):

    # required fields

    subject_id = ReferenceField('subject', required=(True))
    treatment_type = StringField(choices=Treatment_types, required=(True))
    treatment_date = DateField(unique_with=['subject_id', 'treatment_type'], required=(True))

    # optional fields
    outcome_ILAE = ListField(IntField())
    outcome_Engel = ListField(IntField())
    outcome_date = ListField(DateField())

    #optional fields specific to treatment_type = surgery

    surgery_type = StringField()
    surgery_pathology = StringField()
    surgery_hemisphere = StringField(choices=Hemispheres)

## Level Three: Channels

class channel(Document):

    # required fields
    subject_id = ReferenceField('subject', required=(True))
    exam_id = ReferenceField('exam', required=(True))
    channel_name = StringField(required=(True), unique_with=['exam_id'])

    # optional fields

    channel_hemisphere = StringField(choices=Hemispheres)
    electrode_type = StringField(choices=electrode_types)
    location = ListField(FloatField())
    ROI_ids = ListField(IntField())
    ROI_name = ListField(StringField())

    is_resected = BooleanField() # if within 5 mm of resection mask, DP B2
```

```
is_spiking = BooleanField()
is_within_structurally_abnormal_tissue = BooleanField()
is_within_soz = BooleanField()
```

Level Three: icEEG Segments

```
class iceeg_segment(Document):
```

required fields

```
subject_id = ReferenceField('subject', required=(True))
exam_id = ReferenceField('exam', required=(True))
icEEG_datetime = DateTimeField(unique_with=['exam_id'], required=(True))
icEEG_duration = FloatField(required=(True))
icEEG_sampling_frequency = FloatField(required=(True))
icEEG_reference = StringField(choices=icEEG_reference_types, required=(True))
icEEG_wake_status = StringField(choices=sleep_stages)
icEEG_segment_file_location = StringField(required=(True))
```

no optional fields

Level Three: Seizures

```
class seizure(Document):
```

required fields

```
subject_id = ReferenceField('subject', required=(True))
exam_id = ReferenceField('exam', required=(True))
seizure_datetime = DateTimeField(unique_with=['exam_id'], required=(True))
```

optional fields

```
seizure_duration = FloatField()
seizure_type = StringField(choices=seizure_types)
seizure_onset_channels = ListField(StringField())
```

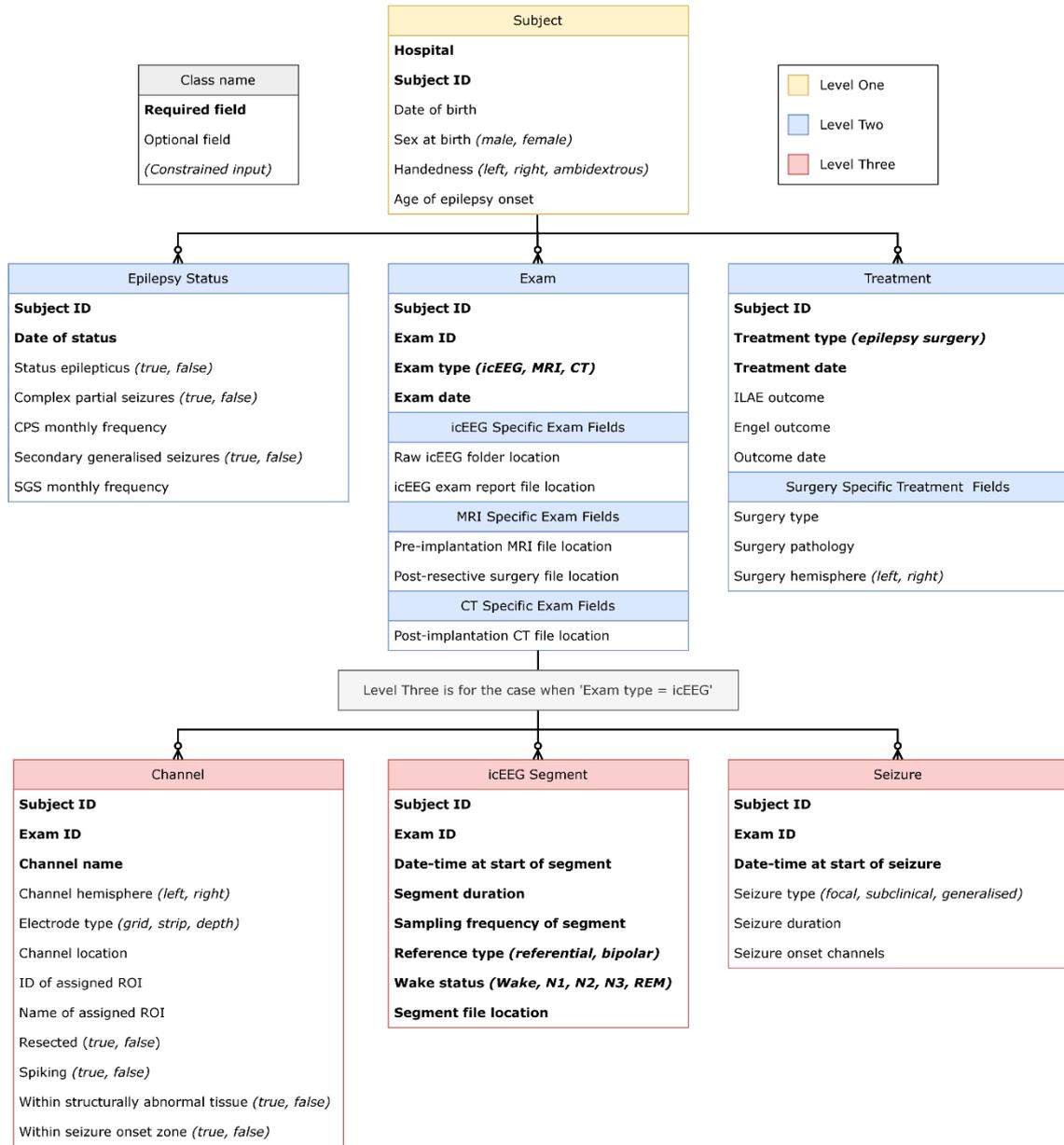


Figure 5: Visualisation of the recommended database structure

D. Normative mapping

The **goal of Pipeline D** is to use the database output to generate a regional map of normative brain activity which encompasses all interictal channel-level data.

Before starting this pipeline, some **Decision Points (DPs)** must be addressed. These are the reader's choice and will affect the resulting normative map. We provide the decisions made in Taylor *et al.* (1) where possible.

- **Decision Point D1:** Will there be a constraint on the type of electrode used in the normative map (i.e., grid, depth)? Certain signal processing properties may differ across electrode types due to their physical properties or location differences. One study finds that signal characteristics appear to be similar, but not identical between electrode types (37). Nevertheless, Taylor *et al.* (1) include all intracranial electrode types in their normative map.
- **Decision Point D2:** Will there be any further precautionary steps to identify noisy channels before constructing the map? If so, what will these be? Taylor *et al.* (1) locate and remove channels with outlier signal range and/or variance relative to the other channels – a more detailed outline of the process can be found in **Supplementary Information 1**.
- **Decision Point D3:** Which preprocessing steps will be applied to the icEEG segments? Preprocessing choices can impact brain signals, including the referencing approach (38,39). Briefly, Taylor *et al.* (1) apply a common average reference within a subject and bandpass filter. They also downsample the data to a common sampling frequency across all subjects, regions and channels and remove power line noise. Specifics can be found in **Supplementary Information 2**.
- **Decision Point D4:** Which metric(s) will be used for the normative map? Some metrics can be calculated directly from the icEEG segment, such as PSD and phase-amplitude coupling. The PSDs can also be used to calculate further metrics, such as relative band power (RBP), absolute band power, and power shift. Taylor *et al.* (1) construct a normative map of log(RBP) in five frequency bands of interest. See **Supplementary Information 3** for their steps in more detail.
- **Decision Point D5:** When constructing a regional map, there may be instances where a subject has two or more channels assigned to a single ROI – how should the data from multiple channels be consolidated to the ROI level? Taylor *et al.* (1) take the mean of the RBP values to give each subject a single value of RBP per frequency band per implanted region.
- **Decision Point D6:** How will outliers be identified at the regional level? Outliers could skew the normative distributions and can impact downstream results. For example, calculated abnormalities (Pipeline E) can be drastically underestimated due to outliers artificially inflating the spread of the normative distribution. In our ongoing projects, we use an iterative Leave-One-Out approach, by assessing each subject against the rest of the normative map, see **Supplementary Information 4** for the specific method.
- **Decision Point D7:** Following all data selection steps, some subjects may have more than one segment remaining – how will multiple segments originating from the same subject be organised? This decision is

dependent on the research question; the motivation for selecting multiple segments in our ongoing projects is to maximise the chance of a subject still having a segment if one (or more) was labelled as an outlier, and to establish temporal variability/stability (5). Hence, following all data refinement, any remaining ‘backup’ segments can be excluded. The reader might want to handle subjects with multiple segments differently if, for example, segments are to be used for validation or for investigating normative map variability.

Once **Decision Points D1-7** have been finalised, begin Normative Mapping (Pipeline D). This pipeline is applied to *all subjects*. Alternatively, a subset can be held out for Abnormality Mapping (Pipeline E). This pipeline assumes the reader has constructed the database outlined in Pipeline C. We recommend creating a new folder to store outputs from this pipeline, see **Error! Reference source not found.**

1. Query the database created in Pipeline C to obtain the subset of subjects which will be used for normative mapping, applying any desired constraint on electrode type (**DP D1**).

Quality control: Remember to consider, and be consistent with, the decisions made in Temporal Processing (Pipeline A, **DP A1-7**).

2. Load the subset of normative mapping subjects identified in the previous step into the reader’s preferred programming language.
3. Refine the data by removing all channels which are thought to represent pathology (non-normative). If available, this information was recorded in the previous pipelines. Non-normative channels are those which were labelled as:
 - a) Resected
 - b) Spiking
 - c) Within structural abnormalities
 - d) Within the SOZ
4. Refine the data further by removing any channels which are noisy or unsuitable. These are channels which:
 - a) Were labelled as ‘unsuitable’ channels during inspection and extraction of icEEG segments
 - b) Did not localise to an ROI in Pipeline B/do not have mapping information available
 - c) Are algorithmically detected as noisy (**DP D2**)
5. Apply icEEG preprocessing steps to all remaining channels (**DP D3**).

Quality control: Power line noise can differ by country. For example, data from UK hospitals has power line noise at 50 Hz, whilst data from the USA has it at 60 Hz. If removing power line noise, make sure to apply the correct filter to each batch of data.

Quality control: If the reader is considering larger frequency ranges, remember to apply a notch filter for harmonics where applicable.

6. Save the preprocessed segment(s) in a suitable folder location. We recommend one folder per subject, see **Error! Reference source not found.**

7. Calculate metrics for each channel in each preprocessed segment (**DP D4**).
8. Save the metrics calculated in the previous step in a suitable folder location. We again recommend one folder per subject, see **Error! Reference source not found.**

Pause point: At this point, the reader should have preprocessed interictal segment(s) for each subject, with metrics calculated for each. These segments should be stored in an easily accessible, organised manner with non-normative and unsuitable channels removed.

9. Map channel data to the regional level, with consideration of whether to/how to attain a single metric per subject per region (**DP D5**).
10. Check for and label outliers within the current normative data table (**DP D6**).
11. Reduce the normative data by removing all outliers.
12. Reduce the normative data, if necessary, to retain the desired number of segments per subject (**DP D7**).

Quality control: Depending on the intended research question, the reader might check the subject sample size per ROI here, or the distribution of samples, and potentially exclude any regions that are not adequate. For example, suppose the intention is to use the normative map as a baseline to calculate z-scores. In that case, each ROI should have at least 30 samples, and the distribution of samples should be approximately Normal.

13. Save the resultant data table in a suitable location.

Endpoint: A data table which provides a regional, normative map of brain activity using the reader's desired metric(s). The map has had samples thought to represent pathology removed and has been algorithmically tested for outliers. As an example endpoint, Taylor *et al.* (1) take the normative distribution of RBP in a region (in a particular frequency band) to be the distribution of RBPs of all subjects with coverage in that region.

Completion of Pipeline D concludes the fundamental sections of our protocol. The normative map is now ready for downstream analysis.

E. Abnormality mapping (optional)

The **goal of this pipeline** is to compare new subjects to the normative map, with the aim of identifying abnormalities that may represent pathology.

Before starting this pipeline, some **Decision Points (DPs)** must be addressed. These are the reader's choice and will affect the resulting abnormality map. We provide the decisions made in Taylor *et al.* (1) as an example.

- **Decision Point E1:** What method will be used to score subjects against the normative map? To estimate regional abnormality compared to the normative distribution, Taylor *et al.* (1) compute the absolute z-scores of

RBP values. They then compute the maximum RBP across frequency bands for each region, thereby creating a subject-specific abnormality map. For more details, see **Supplementary Information 5**. An alternative could be centile mapping, and channel-wise, as well as region-wise, abnormality mapping is also possible.

- **Decision Point E2 (optional, require resection mask):** What threshold will be implemented to define a resected region? For each subject in the abnormality cohort, Taylor *et al.* (1) define a region as resected if >25% of channels within the region were removed, otherwise, regions were considered spared.
- **Decision Point E3 (optional, requires resection mask):** How will the difference in abnormality scores between resected and spared regions be quantified? Abnormality measures can be computed for resected and spared regions, allowing for direct comparison. Taylor *et al.* (1) use the distinguishability statistic (D_{RS}), which is the area under the receiver operating curve. See **Supplementary Information 6** for a more detailed explanation.

Additionally, some of the decisions made in Pipeline D (when constructing the normative map) need to be maintained for Abnormality Mapping (**DP D1-5**), otherwise, the comparison between identified abnormalities and the normative baseline will not be useful. The preprocessing steps (**DP D3**) and chosen metric (**DP D4**) must be the same as in Pipeline D for abnormality mapping to be viable. We recommend keeping the electrode type (**DP D1**) and method for handling multiple channels per subject per region (**DP D5**) the same also. Noisy channel detection (**DP D2**) should be applied here but strictly to identify non-physiological channels, it should not risk excluding any genuine interictal spiking, for example.

Regional outlier detection (**DP D6**) should not be implemented as this may identify the abnormalities that are actively being sought in this pipeline. For simplicity, this pipeline assumes only one segment per subject is retained for abnormality mapping, negating the need for **DP D7**.

Once **Decision Points E1-3** have been finalised, begin Abnormality Mapping (Pipeline E). Abnormality mapping is either done using a subset of subjects held out during Normative Mapping (Pipeline D) or using *new subjects*. Unseen subjects will need to go through Pipeline A-C before this one. This pipeline assumes the reader has the normative map created in Pipeline D ready to use.

1. Query the database created in Pipeline C to obtain the subset of subjects which will be used for abnormality mapping, applying any desired constraint on electrode type (**DP D1**).

Quality control: Remember to consider, and be consistent with, the decisions made in Temporal Processing (Pipeline A, **DP A1-7**).

2. Load the subset of abnormality mapping subjects identified in the previous step into the reader's preferred programming language.
3. Refine the data by removing any channels which are noisy or unsuitable. These are channels which:
 - a) Were labelled as 'unsuitable' channels during inspection and extraction of icEEG segments
 - b) Did not localise to an ROI in Pipeline B/do not have mapping information available
 - c) Are algorithmically detected as noisy (**DP D2**)

Note: We do not remove suspected pathological channels (as in step 3 of Pipeline D), the goal is no longer to obtain normative data

4. Apply icEEG preprocessing steps to all remaining channels (**DP D3**)

Quality control: Power line noise can differ by country. For example, data from UK hospitals has power line noise at 50 Hz, whilst data from the USA has it at 60 Hz. If removing power line noise, make sure to apply the correct filter to each batch of data.

Quality control: If the reader is considering larger frequency ranges, remember to apply a notch filter for harmonics where applicable.

5. Save the preprocessed segment(s) in a suitable folder location.

Note: We recommend one folder per subject, see **Error! Reference source not found.** for the organisation of normative mapping files. Similar organisation can be implemented for abnormality mapping files.

6. Calculate metrics for each channel in each segment (**DP D4**)
7. Save the metrics calculated in the previous in a suitable folder location.

Pause point: At this point, the reader should have preprocessed interictal segment(s) for each subject, with metrics calculated for each. These segments should be stored in an easily accessible, organised manner with unsuitable channels removed. Note, so far we are following the same steps as Pipeline D, but have not excluded non-normative data.

8. Map channel data to the regional level, with consideration of whether to/how to attain a single metric per subject per region (**DP D5**)
9. Quantify regional abnormalities through comparison to the normative map for each subject (**DP E1**)
10. [optional] Identify which regions were resected and which regions were spared for each subject (**DP E2**)
11. [optional] Compute abnormality measures between resected and spared regions (**DP E3**)

End point: A regional abnormality map for each subject in the abnormality cohort, highlighting areas which may be pathological through comparison to the normative map.

Data analysis

This protocol aims to prepare neuroscience researchers for the data analysis stage by equipping them with the tools to construct normative maps. It does not offer advice on the data analysis itself.

However, our lab's published works provide some examples of the subsequent data analysis stage, demonstrating the type of results which can be obtained using this normative mapping protocol (or aspects of it). Below we highlight an

example result from some of our research articles, see the ‘Methods’ section of each reference for details on the data analysis procedures.

Pipelines A-E

- Taylor *et al.* (1) use normative and abnormality maps to demonstrate that interictal icEEG can localise epileptogenic tissue (Figure 3, Figure 4)
- Kozma *et al.* (2) find that relative complete icEEG band power is more effective for distinguishing between surgical outcome groups than the periodic or aperiodic component alone (Figure 4)
- Horsley *et al.* (3) demonstrate that both connectivity and icEEG abnormalities can localise epileptogenic tissue (Figure 4), and that dMRI abnormalities could inform icEEG electrode placement (Figure 5)
- Owen *et al.* (4) find that a combination of magnetoencephalography and icEEG abnormalities is predictive of surgical outcome (Figure 4)
- Wang *et al.* (5) show that for icEEG, the D_{RS} statistic (a measure of how different surgically resected and spared tissue is) is relatively consistent over time and able to separate surgical outcome groups (Figure 2, Figure 3)
- Steinschneider *et al.* (6) use icEEG from a child with left insular ganglioglioma to identify abnormalities by comparison to a normative map (Figure 5, Table 3)

Pipelines A-D

- Janiukstyte *et al.* (7) show that normative maps generated from scalp EEG correlate positively with those generated from icEEG and magnetoencephalography (Figure 3)

Pipeline B

- Wang *et al.* (8) show that normalising for spatial proximity between nearby icEEG electrodes improves the prediction of post-surgery seizure outcomes (Figure 2)
- Thornton *et al.* (9) use long-term icEEG recordings to show that circadian and ultradian rhythms are diminished in pathological brain tissue (Figure 1, Figure 2)
- Gascoigne *et al.* (10) construct resection masks to demonstrate that the incomplete resection of the icEEG SOZ is not associated with surgery outcomes (Figure 2)

Pipeline B – resection masks

- Taylor *et al.* (11) use pre- and post-surgery MRIs to infer that the impact of surgery leads to a reduction of <10% in efficiency in the majority of subjects (Figure 4)
- Owen *et al.* (12) find that markers of surgical failure mechanisms, such as failure to resect magnetoencephalography abnormalities, discriminate between surgical outcome groups (Figure 5)

Validation of protocol

This protocol (or aspects of it) has been validated in the twelve published articles outlined in the section above (**Data Analysis**). The DOIs of the validation articles are also listed after the **Keywords** section.

General notes

Throughout the protocol, **Decision Points** are included where appropriate. These are intended to be methodological choices which impact the resultant normative map and subsequent analysis. They are parameters which, when changed, lead to a different outcome, for example, constructing a normative map of wake activity versus REM activity. Decisions which should not impact outcomes are not included. For instance, whether the reader extracts interictal segments using MATLAB, R or Python, should not impact downstream results.

The inclusion of **Decision Points** allows the reader to tailor the protocol to their own research question. It also facilitates repeatability, allowing the protocol to be implemented for various research questions by merely changing one (or more) **Decision Point(s)**. Nevertheless, as a starting point for most **Decision Points**, we include the choices made in our primary research paper (1). If the reader were to follow the protocol through with those decisions implemented, they would be in a position to reproduce the results in that paper.

Acknowledgements

We thank members of the Computational Neurology, Neuroscience & Psychiatry Lab (www.cnp-lab.com) for discussions on the analysis and manuscript. We also thank each author on our list of published validation papers.

Lastly, we thank our funders. H.W, S.J.G and C.S are funded by the Engineering and Physical Sciences Research Council (EP/L015358/1). Y.W and P.N.T are funded by the UKRI Future Leaders Fellowships (MR/V026569/1, MR/T04294X/1). N.E is funded by Epilepsy Research Institute UK.

Competing interests

The authors have no competing interests to disclose.

Ethical considerations

Any work involving data from human subjects must go through the proper channels and attain ethical approval. This approval should be included in any subsequent research. For example, we have obtained NHS database ethics to ensure full anonymisation and databasing across NHS sites (40) separately to permission to handle and analyse anonymised data from the University Ethics Committee at Newcastle University (ref. 12721/2018). In terms of data anonymisation, care should be taken to ensure subjects with rare, identifiable epilepsy conditions, for example, cannot be identified within a cohort. Formal guidelines on databasing at scale to ensure anonymity exist elsewhere (41–43)

As previously mentioned, public datasets are also available. Some of our published works (1–7) use the RAM cohort (<http://memory.psych.upenn.edu/RAM>).

References

1. Taylor PN, Papasavvas CA, Owen TW, Schroeder GM, Hutchings FE, Chowdhury FA, et al. Normative brain mapping of interictal intracranial EEG to localize epileptogenic tissue. *Brain*. 2022 Apr 29;145(3):939–49.
2. Kozma C, Schroeder G, Owen T, De Tisi J, McEvoy AW, Miserocchi A, et al. Identifying epileptogenic abnormality by decomposing intracranial EEG and MEG power spectra. *Journal of Neuroscience Methods*. 2024 Aug;408:110180.
3. Horsley JJ, Thomas RH, Chowdhury FA, Diehl B, McEvoy AW, Miserocchi A, et al. Complementary structural and functional abnormalities to localise epileptogenic tissue. *eBioMedicine*. 2023 Nov;97:104848.
4. Owen TW, Janiukstyte V, Hall GR, Chowdhury FA, Diehl B, McEvoy A, et al. Interictal magnetoencephalography abnormalities to guide intracranial electrode implantation and predict surgical outcome. *Brain Communications*. 2023 Nov 1;5(6):fcad292.
5. Wang Y, Schroeder GM, Horsley JJ, Panagiotopoulou M, Chowdhury FA, Diehl B, et al. Temporal stability of intracranial electroencephalographic abnormality maps for localizing epileptogenic tissue. *Epilepsia*. 2023;64(8):2070–80.
6. Steinschneider M, Rhone AE, Taylor PN, Nourski KV, Dlouhy BJ, Howard MA. Insights into epileptic aphasia: Intracranial recordings in a child with a left insular ganglioglioma. *Epilepsy & Behavior Reports*. 2024;28:100715.
7. Janiukstyte V, Owen TW, Chaudhary UJ, Diehl B, Lemieux L, Duncan JS, et al. Normative brain mapping using scalp EEG and potential clinical application. *Sci Rep*. 2023 Aug 18;13(1):13442.
8. Wang Y, Sinha N, Schroeder GM, Ramaraju S, McEvoy AW, Miserocchi A, et al. Interictal intracranial electroencephalography for predicting surgical success: The importance of space and time. *Epilepsia*. 2020;61(7):1417–26.
9. Thornton C, Panagiotopoulou M, Chowdhury FA, Diehl B, Duncan JS, Gascoigne SJ, et al. Diminished circadian and ultradian rhythms of human brain activity in pathological tissue in vivo. *Nat Commun*. 2024 Oct 2;15(1):8527.
10. Gascoigne SJ, Evans N, Hall G, Kozma C, Panagiotopoulou M, Schroeder GM, et al. Incomplete resection of the intracranial electroencephalographic seizure onset zone is not associated with postsurgical outcomes. *Epilepsia* [Internet]. 2024 Sep [cited 2025 Jan 9];65(9). Available from: <https://onlinelibrary.wiley.com/doi/10.1111/epi.18061>
11. Taylor PN, Sinha N, Wang Y, Vos SB, de Tisi J, Miserocchi A, et al. The impact of epilepsy surgery on the structural connectome and its relation to outcome. *NeuroImage: Clinical*. 2018;18:202–14.

12. Owen TW, Schroeder GM, Janiukstyte V, Hall GR, McEvoy A, Miserocchi A, et al. MEG abnormalities and mechanisms of surgical failure in neocortical epilepsy. *Epilepsia*. 2023;64(3):692–704.
13. Frauscher B, von Ellenrieder N, Zemann R, Doležalová I, Minotti L, Olivier A, et al. Atlas of the normal intracranial electroencephalogram: Neurophysiological awake activity in different cortical areas. *Brain : a journal of neurology*. 2018 Mar;141.
14. Bernabei JM, Sinha N, Arnold TC, Conrad E, Ong I, Pattnaik AR, et al. Normative intracranial EEG maps epileptogenic tissues in focal epilepsy. *Brain*. 2022 May;145(6):1949–61.
15. Groppe DM, Bickel S, Keller CJ, Jain SK, Hwang ST, Harden C, et al. Dominant frequencies of resting human brain activity as measured by the electrocorticogram. *NeuroImage*. 2013;79:223–33.
16. Owen LLW, Muntianu TA, Heusser AC, Daly PM, Scangos KW, Manning JR. A Gaussian Process Model of Human Electrocorticographic Data. *Cerebral Cortex*. 2020 Sep 3;30(10):5333–45.
17. Betzel RF, Medaglia JD, Kahn AE, Soffer J, Schonhaut DR, Bassett DS. Structural, geometric and genetic factors predict interregional brain connectivity patterns probed by electrocorticography. *Nat Biomed Eng*. 2019 May 27;3(11):902–16.
18. Kalamangalam GP, Long S, Chelaru MI. A neurophysiological brain map: Spectral parameterization of the human intracranial electroencephalogram. *Clinical Neurophysiology*. 2020 Mar;131(3):665–75.
19. Frauscher B, Von Ellenrieder N, Zemann R, Rogers C, Nguyen DK, Kahane P, et al. High-Frequency Oscillations in the Normal Human Brain. *Annals of Neurology*. 2018 Sep;84(3):374–85.
20. Frauscher B, Von Ellenrieder N, Dolezalova I, Bouhadoun S, Gotman J, Peter-Derex L. Rapid Eye Movement Sleep Sawtooth Waves Are Associated with Widespread Cortical Activations. *J Neurosci*. 2020 Nov 11;40(46):8900–12.
21. Simpson C, Hall G, Duncan JS, Wang Y, Taylor PN. Automated generation of epilepsy surgery resection masks; The RAMPS pipeline.
22. Arnold TC, Muthukrishnan R, Pattnaik AR, Sinha N, Gibson A, Gonzalez H, et al. Deep learning-based automated segmentation of resection cavities on postsurgical epilepsy MRI. *NeuroImage: Clinical*. 2022;36:103154.
23. Pérez-García F, Dorent R, Rizzi M, Cardinale F, Frazzini V, Navarro V, et al. A self-supervised learning strategy for postoperative brain cavity segmentation simulating resections. *Int J CARS*. 2021 Oct;16(10):1653–61.
24. Casseb RF, De Campos BM, Morita-Sherman M, Morsi A, Kondylis E, Bingaman WE, et al. ResectVol: A tool to automatically segment and characterize lacunas in brain images. *Epilepsia Open*. 2021 Dec;6(4):720–6.

25. Cahill V, Sinclair B, Malpas CB, McIntosh AM, Chen Z, Vivash LE, et al. Metabolic patterns and seizure outcomes following anterior temporal lobectomy. *Annals of Neurology*. 2019 Feb;85(2):241–50.
26. Courtney MR, Sinclair B, Neal A, Nicolo JP, Kwan P, Law M, et al. Automated segmentation of epilepsy surgical resection cavities: Comparison of four methods to manual segmentation. *NeuroImage*. 2024 Aug;296:120682.
27. Fischl B. FreeSurfer. *NeuroImage*. 2012 Aug;62(2):774–81.
28. Hamilton LS, Chang DL, Lee MB, Chang EF. Semi-automated Anatomical Labeling and Inter-subject Warping of High-Density Intracranial Recording Electrodes in Electrocorticography. *Front Neuroinform*. 2017 Oct 31;11:62.
29. Tustison NJ, Cook PA, Holbrook AJ, Johnson HJ, Muschelli J, Devenyi GA, et al. The ANTsX ecosystem for quantitative biological and medical imaging. *Sci Rep*. 2021 Apr 27;11(1):9068.
30. Combrisson E, Vallat R, Eichenlaub JB, O'Reilly C, Lajnef T, Guillot A, et al. Sleep: An Open-Source Python Software for Visualization, Analysis, and Staging of Sleep Data. *Front Neuroinform*. 2017 Sep 21;11:60.
31. Von Ellenrieder N, Peter-Derex L, Gotman J, Frauscher B. SleepSEEG: automatic sleep scoring using intracranial EEG recordings only. *J Neural Eng*. 2022 Apr 1;19(2):026057.
32. Lodder SS, Askamp J, Van Putten MJAM. Inter-ictal spike detection using a database of smart templates. *Clinical Neurophysiology*. 2013 Dec;124(12):2328–35.
33. Gabor AJ, Seyal M. Automated interictal EEG spike detection using artificial neural networks. *Electroencephalography and Clinical Neurophysiology*. 1992 Nov;83(5):271–80.
34. Desikan RS, Ségonne F, Fischl B, Quinn BT, Dickerson BC, Blacker D, et al. An automated labeling system for subdividing the human cerebral cortex on MRI scans into gyral based regions of interest. *NeuroImage*. 2006 Jul;31(3):968–80.
35. McCarthy CS, Ramprashad A, Thompson C, Botti JA, Coman IL, Kates WR. A comparison of FreeSurfer-generated data with and without manual intervention. *Front Neurosci* [Internet]. 2015 Oct 21 [cited 2025 Feb 3];9. Available from: <http://journal.frontiersin.org/Article/10.3389/fnins.2015.00379/abstract>
36. FsTutorial/AnatomicalROI - Free Surfer Wiki [Internet]. [cited 2025 Feb 3]. Available from: <https://surfer.nmr.mgh.harvard.edu/fswiki/FsTutorial/AnatomicalROI>
37. Alkawadri CI, Yan Q, Kocuglu Kinal AG, Spencer DD, Alkawadri R. Comparison of EEG Signal Characteristics of Subdural and Depth Electrodes. *Journal of Clinical Neurophysiology* [Internet]. 2024 Dec 30 [cited 2025 Feb 3]; Available from: <https://journals.lww.com/10.1097/WNP.0000000000001139>

38. Michelmann S, Treder MS, Griffiths B, Kerrén C, Roux F, Wimber M, et al. Data-driven re-referencing of intracranial EEG based on independent component analysis (ICA). *Journal of Neuroscience Methods*. 2018 Sep;307:125–37.
39. Shirhatti V, Borthakur A, Ray S. Effect of Reference Scheme on Power and Phase of the Local Field Potential. *Neural Computation*. 2016 May 1;28(5):882–913.
40. The National Epilepsy and Neurology Database [Internet]. [cited 2024 Feb 7]. Available from: <https://www.hra.nhs.uk/planning-and-improving-research/application-summaries/research-summaries/national-epilepsy-neurology-database/>
41. Olatunji IE, Rauch J, Katzensteiner M, Khosla M. A Review of Anonymization for Healthcare Data. *Big Data*. 2024 Dec 1;12(6):538–55.
42. El Emam K, Arbuckle L. *Anonymizing health data: case studies and methods to get you started*. O'Reilly Media, Inc.; 2013.
43. Rodriguez A, Tuck C, Dozier MF, Lewis SC, Eldridge S, Jackson T, et al. Current recommendations/practices for anonymising data from clinical trials in order to make it available for sharing: A scoping review. *Clinical Trials*. 2022 Aug;19(4):452–63.

Supplementary information

1. Algorithmic detection of noisy channels

Below we provide an example method for detecting noisy icEEG channels, which we have used in some of our published works. Its purpose is to find channels with outlier signal range and/or variance relative to the other channels.

Two rounds of detection are performed. The first round is before preprocessing with (by default) less stringent detection thresholds; the second round is after basic preprocessing. All the specified thresholds and preprocessing steps are the reader's choice, here we simply outline our own choices.

1. [Round one] If a channel has a signal variance and/or range greater than 16, mark it as an outlier
2. Apply preprocessing steps:
 - a) Common average reference

- b) A 4th order, zero-phase shift Butterworth filter between 0 and 100 Hz
- c) A notch filter to eliminate the location-specific power line noise (i.e., a 2 Hz window between 49 and 51 Hz to eliminate 50 Hz noise for UK hospitals)

Quality control: If the reader is considering larger frequency ranges, remember to apply a notch filter for harmonics

- 3. [Round two] If a channel has a signal variance and/or range greater than 12, mark it as an outlier

2. icEEG preprocessing steps

icEEG preprocessing steps are the reader's choice. Below we outline the specifics of the preprocessing steps used in Taylor *et al.* (1):

- a) Common average reference
- b) Bandpass filtering between 0.5-80 Hz using a 4th order, zero phase-shift Butterworth filter
- c) Downsampling to 200 Hz with an anti-alias filter
- d) Removal of power line noise using a notch filter with a 2 Hz window (e.g. 49-51 Hz for UK data, 59-61 Hz for USA data)

3. Calculating a normative metric: log(RBP)

The normative mapping metric(s) is reader's choice. Taylor *et al.* (1) construct a map of log(RBP) in five frequency bands of interest. The steps for calculating this metric are:

- 1. Calculate PSDs using Welch's method with a 2-second window and a 1-second overlap
- 2. Calculate the average band power in five standard frequency bands, δ (1-4 Hz), θ (4-8 Hz), α (8-13Hz), β (13-30 Hz) and γ (30-77.5 Hz)

Note: Typically, metrics such as RBP would not be calculated in any frequency ranges which have been removed using a notch filter.

Quality control: In the γ range, Taylor *et al.* (1) exclude the frequency ranges 47.5-52.5 and 57.5-62.5. This is because the data originated from both UK and USA hospitals, and therefore different subjects had different notch filters for power line noise (50 and 60 Hz for UK and USA subjects respectively). Removing those frequencies for all data ensures the map, and any results are not impacted by the differing preprocessing steps.

- 3. Apply a log10 transformation to the band power estimates
- 4. Normalise the band power estimates to sum to one in each channel (L1-normalisation)

4. Identifying outliers at the regional level

Outlier segments can be identified at the regional level to ensure that the final regional map is as normative as possible. Below we outline one possible method for outlier detection:

1. Hold out one regional-level icEEG segment from the remaining normative data
2. Compute the z-score of this segment using all remaining segments as a normative baseline

Quality control: The reader should check the distribution of values (e.g., RBP) follows a Normal distribution when using this method

3. Label the segment as an outlier if $z > 2$ (the threshold is the reader's choice)
4. Repeat steps 1-3 for all segments
5. Hold out all segments labelled as outliers
6. Repeat steps 1-4 up to, for example, ten times, or stop if a round does not identify any new outliers

5. Scoring new subjects against the normative map

To estimate the abnormality of a region's RBP for a new subject, Taylor *et al.* (1) score against the normative map of RBP values by computing the absolute z-score. This is done using the following equation:

$$|z_{i,j}| = \left| \frac{x_{i,j} - \mu_{i,j}}{\sigma_{i,j}} \right|$$

where i represents the ROI, j is the frequency band of interest, x is the RBP value for the subject and μ and σ are the mean and standard deviations of the RBPs in the normative map respectively. The final abnormality value (for a given ROI) is then the maximum absolute z-score across frequency bands. Calculating these values for each of the subject's implanted regions forms the abnormality map, with higher absolute z-scores indicating greater abnormality relative to the normative map. This can be visualised for ease of interpretation – see Taylor *et al.* Figure 2 (1).

6. Comparing resected and spared regions

If resection masks have been created for the abnormality cohort, abnormality measures between resected and spared regions can also be calculated. Taylor *et al.* (1) use the distinguishability statistic (D_{RS}), which is the area under the receiver operating curve, and equivalent to the normalized Mann–Whitney U-statistic (1,8). A D_{RS} value > 0.5 indicates that spared regions were more abnormal (higher absolute z-score) than resected regions, whereas D_{RS} values < 0.5 indicates the opposite—i.e. resected regions were more abnormal. See Taylor *et al.* Figure 3 for a visualisation (1).