

# M<sup>3</sup>Builder: A Multi-Agent System for Automated Machine Learning in Medical Imaging

Jinghao Feng<sup>\*,1,2</sup> Qiaoyu Zheng<sup>\*,1,2</sup> Chaoyi Wu<sup>1,2</sup> Ziheng Zhao<sup>1,2</sup>  
Ya Zhang<sup>1,2</sup> Yanfeng Wang<sup>†,1,2</sup> and Weidi Xie<sup>†,1,2</sup>

<sup>1</sup>Shanghai Jiao Tong University

<sup>2</sup>Shanghai AI Laboratory

\* Equal Contribution


† Corresponding Author

**Abstract.** Agentic AI systems have gained significant attention for their ability to autonomously perform complex tasks. However, their reliance on well-prepared tools limits their applicability in the medical domain, which requires to train specialized models. In this paper, we make three contributions: **(i)** We present **M<sup>3</sup>Builder**, a novel multi-agent system designed to automate machine learning (ML) in medical imaging. At its core, M<sup>3</sup>Builder employs four specialized agents that collaborate to tackle complex, multi-step medical ML workflows, from automated data processing and environment configuration to self-contained auto debugging and model training. These agents operate within a medical imaging ML workspace, a structured environment designed to provide agents with free-text descriptions of datasets, training codes, and interaction tools, enabling seamless communication and task execution. **(ii)** To evaluate progress in automated medical imaging ML, we propose **M<sup>3</sup>Bench**, a benchmark comprising four general tasks on 14 training datasets, across five anatomies and three imaging modalities, covering both 2D and 3D data. **(iii)** We experiment with seven state-of-the-art large language models serving as agent cores for our system, such as Claude series, GPT-4o, and DeepSeek-V3. Compared to existing ML agentic designs, M<sup>3</sup>Builder shows superior performance on completing ML tasks in medical imaging, **achieving a 94.29% success rate using Claude-3.7-Sonnet as the agent core**, showing huge potential towards fully automated machine learning in medical imaging.

**Keywords:** Agentic System · Medical Imaging · Machine Learning.

## 1 Introduction

Large Language Model (LLM)-powered agentic systems have demonstrated remarkable success across diverse domains, leveraging their ability to orchestrate specialized tools and solve complex, multi-step tasks with precision. However, their application in medical domain remains challenging, due to a shortage of well-prepared tools. Such challenge arises from two aspects: first, the complexity of medical workflows—spanning diverse diseases, imaging modalities, and task requirements—makes tool development and integration difficult; second, while

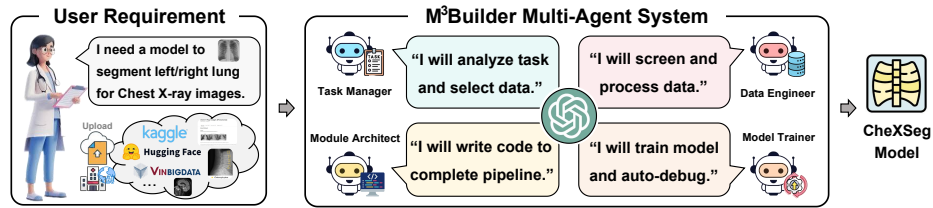

**Fig. 1.** General Pipeline of M<sup>3</sup>Builder. Based on user requirements and available raw clinical datasets, four agents play different roles and collaborate to complete the procedure from task analysis, data processing to training model.

clinicians are experts in their field, they often lack the coding skills or resources to create tools, widening the gap between clinical needs and tailored AI solutions.

In this paper, we propose **M<sup>3</sup>Builder** (M<sup>3</sup> stands for **M**ulti-agent, **M**achine Learning, and **M**edical Imaging), a novel agentic framework for automating machine learning (ML) workflows in medical imaging. At its core, M<sup>3</sup>Builder employs a multi-agent collaboration framework, where four specialized role-playing LLM agents work together to execute complex, multi-step workflows, collaboratively managing task requirement understanding, procedural planning, raw data processing, environment configuration, automated self-correction and model training execution.

Supporting this framework is the Medical Imaging ML workspace, a structured environment that provides the agents with initial datasets, pipeline templates, and interaction tools, enabling seamless human-in-the-loop. Together, the agents and the workspace form a cohesive system capable of autonomously developing AI tools for medical imaging analysis, bridging gaps in existing toolsets and facilitating the creation of a dynamic agentic ecosystem capable of evolving with minimal human intervention. Given a medical imaging-specific ML task and raw training data, **M<sup>3</sup>Builder** autonomously manages the entire development process, from data preparation to model construction, training, and deployment.

In addition, we introduce a new benchmark, **M<sup>3</sup>Bench**, designed to evaluate the performance for automated machine learning in medical imaging. **M<sup>3</sup>Bench** comprises four general tasks: *organ segmentation*, *anomaly detection*, *disease diagnosis*, and *report generation*. These tasks are further associated with 14 detailed training datasets, spanning five anatomical regions and three primary imaging modalities, encompassing both 2D and 3D models. By covering a broad range of tasks and datasets, M<sup>3</sup>Bench provides a standardized foundation for assessing the capabilities of agentic systems in automated ML in medical imaging.

Experimentally, we compare seven state-of-the-art (SOTA) LLMs as the agent core for M<sup>3</sup>Builder, including GPT-4o [10], Claude-Sonnet series [1], and DeepSeek-V3 [17]. Later, to validate the effectiveness of the proposed workspace and multi-agent collaboration framework, we compare against other state-of-the-art ML agentic systems, including ML-AgentBench [9], Aider [6], Cursor Composer [24], Windsurf Cascade [3] and Copilot Edits [7]. Our results demonstrate

that **M<sup>3</sup>Builder** consistently achieves superior performance across a range of metrics, advancing the state-of-the-art in fully autonomous ML workflows for medical imaging, **achieving a 94.29% success rate using Claude-3.7-Sonnet as the agent core**, with 86.67% on organ segmentation, 100% on anomaly detection, 95% on disease diagnosis, and 93.33% on report generation.

## 2 Method

In this section, we present the proposed **M<sup>3</sup>Builder** in detail, starting with the problem formulation, then a detailed description on the ML workspace initialization and the multi-agent collaboration framework.

### 2.1 Problem Formulation

Given a task description on medical imaging analysis, denoted as  $\mathcal{T}$ , our objective is to automatically construct a functional AI model via multi-agents collaboration. As shown in Fig. 2, our proposed framework **M<sup>3</sup>Builder** comprises two key components: the medical imaging ML workspace ( $\mathcal{W}$ ), and the multi-agent collaboration framework ( $\mathcal{A}$ ). Specifically, the workspace includes three elements: *data cards*, *toolset descriptions*, and *code templates*. The data cards are represented in natural language, while the toolset descriptions and code templates are provided in `Python`. Together, these elements form a structured ‘environment’ that guides the automatic AI workflow.

Building on the workspace ( $\mathcal{W}$ ), the multi-agent framework composes of four LLM agents with distinct roles, *i.e.*,  $\mathcal{A} = \{a_1, a_2, a_3, a_4\}$ . These agents adopt a divide-and-conquer strategy to collaboratively address the complexities of the AI task. The framework iteratively performs code generation, editing and executing using tools defined by *toolset descriptions*, until a functional AI model is successfully produced. This process can be expressed as:


$$\{\mathcal{C}_i, \mathcal{R}_i\} = \mathcal{A}(\mathcal{C}_{i-1}, \mathcal{R}_{i-1}, \mathcal{T} \mid \mathcal{W}), \quad (1)$$

where  $\mathcal{C}_i$  denotes code and scripts generated or edited in the  $i_{th}$  iteration,  $\mathcal{R}_i$  denotes the compiler feedback in `Python` environment, with  $\mathcal{C}_0 = \mathcal{R}_0 = \emptyset$ .

In the following sections, we present more details for the AI workspace initialization and the multi-agent collaboration framework.

### 2.2 Workspace Initialization

The ML Workspace for medical imaging analysis is designed to serve three key purposes: (i) it provides the multi-agent framework with metadata of the available datasets, enabling informed dataset selection; (ii) it supplies initial code templates, offering a standardized starting point for agents and demonstrating a typical AI model training pipeline; (iii) it defines and describes all tools available to the agents, restricting their action space to a predefined, complete set of

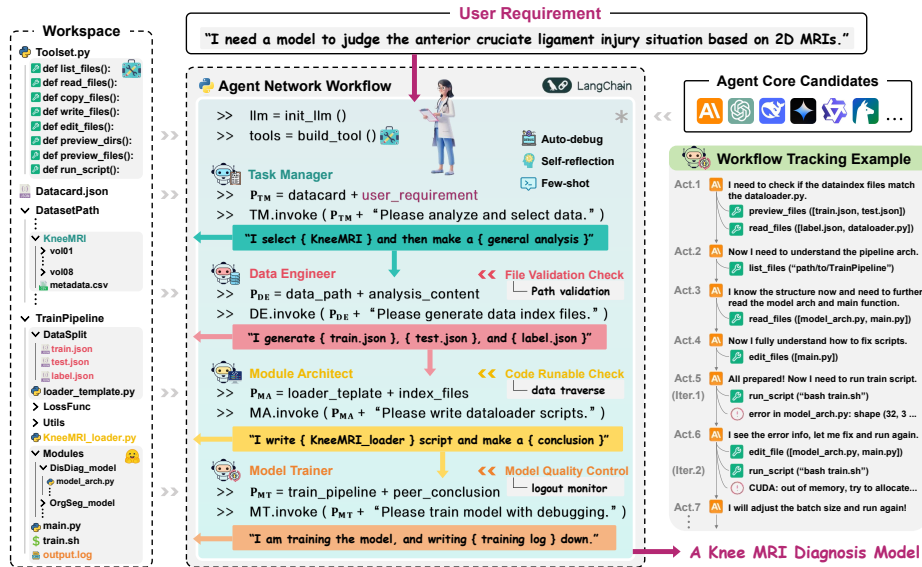

**Fig. 2.** Visualization of M<sup>3</sup>Builder workflow: From user’s free-text request to model delivery. The system integrates user requirements, a workspace with candidate data, tools, and code templates, a network of four specialized collaborative agents. A sample log tracks the Model Trainer agent’s activities during diagnosis model development.

operations. In summary, the workspace equips the multi-agent framework with necessary resources for dataset preparation, a coding foundation, and a well-defined action space, corresponding to dataset descriptions, code templates, and toolset definitions.

**Dataset Preparation.** The workspace includes a range of medical imaging datasets, each accompanied by a structured datacard for standardized descriptions. Each datacard contains key information such as the dataset name, a concise summary (covering data sources, categorical classifications, and scope), and detailed metadata. This flexible format ensures that any description meeting these criteria qualifies as a valid datacard, including documentation provided with the dataset itself. To enable seamless integration, users can add new datasets by completing corresponding datacards. Initially, the workspace provides 14 medical imaging datasets with their datacards as examples.

**Code Template Design.** To streamline AI model training while maintaining flexibility, we prepare a set of standardized code templates based on the Transformers Trainer framework and nnU-Net[11]. These templates are tailored to four primary medical imaging tasks: *disease diagnosis*, *organ segmentation*, *anomaly detection*, and *report generation*. Each task is implemented as a modular package with configurable components, such as the main forward architecture and network backbone options (2D/3D models). Shared features across tasks include a unified selection of loss functions, data augmentation strategies, training util-

ities, and architectural frameworks. By offering these templates, we reduce the complexity of free-form coding while preserving adaptability for diverse AI tasks.

**Interaction Toolset.** To facilitate the agentic system to interact with the PC environment, we develop a comprehensive interaction toolset, comprising eight functions, as shown in the top-right corner of Figure 2. Inspired by MLAgentBench, the toolset includes: *list\_files*, *read\_files*, *copy\_files*, *write\_files*, *edit\_files*, and *run\_script* functions. In addition, we also introduce two specialized tools: *preview\_dirs* and *preview\_files*. The former processes large dataset directories containing numerous files that exceed the LLM’s context window, while the latter extracts key segments from oversized metadata files that cannot be fully processed by the *read\_files* function.

### 2.3 Multi-Agent Collaboration Framework

This section introduces our multi-agent collaboration framework ( $\mathcal{A}$ ), which decomposes the AI task into four sub-tasks and assigns them to four specialized role-playing LLMs agents: Task Manager, Data Engineer, Module Architect, and Model Trainer, denoted as  $\{a_1, a_2, a_3, a_4\}$ , respectively, as illustrated in Fig. 2. Each agent is responsible for a specific role, and together they collaborate iteratively to construct the final AI model. We utilize a set of system prompts to define the role and working logic of each agent, shown in supplementary files.

**Task Manager** acts as the coordinator of the framework. Its primary responsibilities include selecting the most suitable dataset for the task, or alternatively, asking users to upload raw datasets with associated datacard to the workspace as the supplementary of pre-existing datasets, and generating a comprehensive planning document  $\mathcal{P}$  to guide the collaboration among the other agents. Specifically, given a user-provided task description, as exemplified by “user requirements” in Fig. 2, the Task Manager will identify and select the optimal dataset ( $\mathcal{D}$ ) for model training and generate the planning documents. This process can be formally represented as:

$$\{\mathcal{D}, \mathcal{P}\} = a_1(\mathcal{T} \mid \mathcal{W}_d), \quad (2)$$

where  $\mathcal{W}_d$  denotes the data card in the pre-defined workspace.

**Data Engineer** is responsible for dataset preparation and processing. It transforms raw data into a format suitable for model training by performing tasks such as pre-screening the organizational structure of large-scale datasets, analyzing metadata files to extract relevant information, and splitting datasets into training and testing subsets. A critical aspect of the Data Engineer’s role is its iterative interaction with the external compiler environment. It generates, edits, and refines code, incorporating feedback from the compiler, until the code executes successfully. This iterative process ensures the dataset preparation code is both robust and functional. The process can be expressed as:

$$\{\mathcal{C}_i, \mathcal{R}_i\} = a_2(\mathcal{C}_{i-1}, \mathcal{R}_{i-1}, \mathcal{T} \mid \mathcal{D}, \mathcal{P}), \quad (3)$$

where  $\mathcal{C}_i$  represents the  $i^{\text{th}}$  version of the code generated by the Data Engineer, and  $\mathcal{R}_i$  denotes the compiler feedback from the Python environment. Similarly as defined in **Problem Formulation**,  $\mathcal{C}_0 = \mathcal{R}_0 = \emptyset$ .

**Module Architect.** Upon the dataset preparation, the Module Architect integrates essential components into the training pipeline, including developing dataloader scripts, designing appropriate model architecture, initializing the entrance function for training, and selecting other components, such as loss functions, data augmentation strategies, and training utilities. Notably, the Module Architect will iteratively validate the dataloader to ensure it outputs batches with correct shapes and formats. The process can be formulated as:

$$\{\mathcal{C}_i, \mathcal{R}_i\} = a_3(\mathcal{C}_{i-1}, \mathcal{R}_{i-1}, \mathcal{T} \mid \mathcal{D}, \mathcal{W}_c), \quad (4)$$

where  $\mathcal{W}_c$  represents the code templates in the workspace. Similar to the Data Engineer, the code generation process in the Module Architect is also iterative and we denote the final iteration output as  $\mathcal{C}_{\text{MA}}$ . After the integration of these modules, the architect finally synthesizes a summary  $\mathcal{S}$  of all completed work.

**Model Trainer** finalizes the debugging and optimizing the training procedure. Building upon the pipeline established by the Module Architect, the Model Trainer first verifies the completeness and correctness of the training framework thus far. It then selects hyperparameters to meet the model’s specific training requirements, while retaining the authority to modify any part of the code—including model code, dataloader code, and training scripts—based on errors encountered during training. The entire workflow can be expressed as:

$$\{\mathcal{C}_i, \mathcal{R}_i\} = a_4(\mathcal{C}_{i-1}, \mathcal{R}_{i-1}, \mathcal{T} \mid \mathcal{D}, \mathcal{C}_{\text{MA}}, \mathcal{S}). \quad (5)$$

After iteratively performing the above code generation pipeline until successfully executed, the final code will produce a desired AI model.

### 3 Benchmark & Experiments

We evaluate **M<sup>3</sup>Builder**, using 14 diverse ML model development tasks across 4 radiology domains, paired with matched clinical datasets to enable comprehensive testing. The evaluation employs 7 leading LLMs as agent cores: GPT-4 [10], Claude-3.7-Sonnet [1], Claude-3.5-Sonnet [1], DeepSeek-v3 [17], Gemini-2.0-flash [25], Qwen-2.5-max [26], and Llama-3.3-70B [5]. Quantitatively, we assess the system’s effectiveness through analysis on task completion, framework superiority, and agent role-specification.

#### 3.1 Task Definition & Data Preparation

In this paper, we experiment with 14 tasks spanning *organ segmentation*, *anomaly detection*, *disease diagnosis*, and *report generation*. These tasks are systematically categorized by anatomic regions (head & neck, chest, abdomen & pelvis,

**Table 1.** Task Completion Performance Across LLMs. Each experiment undergoes multi-runs, with results shown as successful completions over total rounds (a/b format). **Green** cells indicate that all runs passed, **Yellow** indicates partially passed, and **Red** indicates that all runs failed.

Task	Dataset	LLMs						
		Son.3.7	Son.3.5	GPT4o	QwMax	DeepSeekV3	Gemini2	Llama3.3
OrgSeg	BTCV	5/5	4/5	5/5	3/5	2/5	0/5	1/5
	Verse	4/5	4/5	5/5	2/5	0/5	0/5	0/5
	OASIS	4/5	3/5	3/5	1/5	1/5	0/5	0/5
AnoDet	COV19	5/5	5/5	4/5	3/5	1/5	1/5	0/5
	INS22	5/5	5/5	4/5	2/5	2/5	0/5	0/5
	Panc	5/5	5/5	5/5	4/5	2/5	0/5	1/5
	XDet10	5/5	5/5	5/5	3/5	1/5	1/5	0/5
DisDiag	ADNI	5/5	4/5	2/5	1/5	1/5	0/5	0/5
	KneeMR	4/5	5/5	4/5	2/5	1/5	0/5	0/5
	CC-CCII	5/5	5/5	5/5	0/5	0/5	0/5	0/5
	KidCT	5/5	5/5	4/5	2/5	1/5	0/5	0/5
RepGene	CT-RATE	5/5	4/5	4/5	2/5	1/5	0/5	0/5
	GenomBra	5/5	5/5	4/5	3/5	1/5	0/5	1/5
	IU-Xray	4/5	4/5	3/5	1/5	2/5	0/5	0/5
<b>Average(%)</b>		<b>94.29</b>	90.00	81.43	41.43	22.86	4.29	4.29

limb, spine), imaging modalities (X-ray, CT, MRI), and dimensionality (2D/3D). Each task is precisely defined, for instance: “please build a model for covid-19/pneumonia classification from 3D chest CT images.”

To facilitate model development, we have curated 14 domain-specific datasets: ADNI [12], KneeMRI [23], CC-CCII [27], CT-Kidney [28], BTCV [13], MSD Pancreas [2], VerSe [22,19,16], L2R-OASIS [21], COVID-19 [20], CT-RATE [8], INSTANCE2022 [15], ChestX-Det10 [18], RadGenome-Brain-MRI [14], and IU-Xray [4]. Each dataset is accompanied by a datacard with metadata from their original publications. Additionally, we expand the datacard pool through GPT-generated synthetic entries to simulate real-world scenarios where suitable datasets must be selected from large-scale repositories for model training.

### 3.2 Task Completion Analysis

We design a specific model-building task for each dataset based on its unique characteristics, resulting in a total of 14 tasks. For each task, seven large language models (LLMs) are run independently five times, with a maximum of 100 actions (tool invocations) allowed per execution. The task completion is defined as the successful training of a model with performance on the test set falling within an acceptable range. As the results shown in Tab.1, the performance of different LLMs exhibit significant variation, with Claude-3.7-Sonnet achieving the highest completion rate of 94.29%, while Gemini2.0 and Llama3.3 only reach 4.29%.

**Table 2.** Framework Comparison with SOTAs and Ablations on System Design using Sonnet. Results are averaged over two runs per task in dataset-level. “w/o Colab” represents single-agent execution, and “Iters” means the self-correction rounds.

Agentic System	Completion Runs (Total) ↑					Average Actions (Iters) ↓			
	Seg.	Det.	Diag.	Gen.	Avg(%)	Seg.	Det.	Diag.	Gen.
MLA-Bench	0(6)	0(8)	0(8)	0(6)	0.00	-(-)	-(-)	-(-)	-(-)
Aider	2(6)	3(8)	2(8)	3(6)	35.71	58.0(6.5)	51.3(7.0)	44.5(4.5)	56.7(5.7)
Cursor Comp	1(6)	3(8)	3(8)	2(6)	32.14	42.0(7.0)	35.7(4.3)	36.3(4.0)	51.5(5.5)
Wsurf Casc	1(6)	3(8)	4(8)	2(6)	35.71	39.0(5.0)	36.0(4.3)	35.3(4.8)	48.5(5.5)
Copilot Edits	2(6)	4(8)	3(8)	2(6)	39.29	48.0(3.5)	45.3(3.5)	44.7(4.0)	46.5(4.5)
<b>Ours</b>	4(6)	7(8)	8(8)	4(6)	82.14	34.5(1.8)	25.29(2.4)	35.0(1.9)	32.0(2.3)
w/o Colab	3(6)	3(8)	3(8)	2(6)	39.29	33.6(4.7)	37.5(4.8)	33.3(4.3)	33.5(3.5)
w/o Debug	2(6)	3(8)	1(8)	0(6)	21.43	30.5(1.0)	24.0(1.0)	33.0(1.0)	-(-)
w/o Reflect	4(6)	7(8)	7(8)	4(6)	78.57	28.3(2.3)	24.9(4.1)	35.4(3.9)	41.3(5.8)
w/o Fewshot	3(6)	4(8)	6(8)	3(6)	57.14	37.3(4.0)	23.5(4.3)	32.3(4.1)	34.0(3.0)

The experimental results demonstrate our M<sup>3</sup>Builder possesses impressive and robust capabilities for medical imaging model training automation, particularly when leveraging LLMs that excel in tool utilization and code generation, which significantly enhances the system’s overall performance.

### 3.3 Compare to State-of-the-art Agentic Systems

We compare M<sup>3</sup>Builder with other agentic systems including MAgent-Bench, Aider, Cursor Composer, Windsurf Cascade, and Copilot Edits (all using Sonnet as the agent core). Each system performed each task twice on our workspace under their built-in framework. As shown in Tab. 2, across radiology tasks (Organ Segmentation, Anomaly Detection, Disease Diagnosis, and Report Generation), MAgent-Bench performed poorly due to insufficient data structure understanding capabilities. Other frameworks achieved only moderate success rates (39.29% max) due to single-agent limitations, required human-in-the-loop confirmation, increasing operational complexity and max iteration constraints. Our M<sup>3</sup>Builder demonstrated superior performance with a 42.85% higher average success rate while requiring fewer action steps and execution iterations.

### 3.4 Ablation Study

Here we present ablation studies on our system design, examining the impact of: single-agent versus multi-agent collaboration, auto-debugging capability, self-reflection mechanisms, and workflow few-shot examples. Results in Tab. 2 indicate that self-reflection has minimal influence on system performance, while auto-debugging proves crucial for successful training. Multi-agent collaboration and well-crafted example instructions also significantly impact performance, with their absence resulting in 42.85% and 25.00% performance gaps, respectively.



**Table 3.** Role-specific agent performance on tasks using Sonnet. “Run”, “Act”, “Iter” and “Tkn” respectively denote “execution rounds”, “actions”, “iterations” and “tokens”.

Task	Task Manager				Data Engineer				Module Architect				Model Trainer			
	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn
Seg.	6/6	1.3	1.0	4.3k	5/6	10.2	1.3	72k	5/6	10.5	2.3	74k	4/6	10.7	3.3	115k
Det.	8/8	2.0	1.0	4.8k	8/8	10.0	1.4	92k	7/8	10.6	2.1	61k	7/8	9.3	2.1	67k
Diag.	8/8	2.0	1.0	4.4k	8/8	8.9	1.4	116k	7/8	11.3	2.0	84k	8/8	9.3	2.0	198k
Gen.	6/6	2.0	1.0	4.2k	6/6	8.7	1.3	66k	5/6	10.8	2.5	91k	5/6	11.2	2.2	70k

### 3.5 Analysis on Different Agent Roles

We evaluate the performance of each role-specific agent in **M<sup>3</sup>Builder**, using distinct success criteria: The Task Managers is assessed on its ability to select appropriate training datasets, the Data Engineer on generating valid data index files with correct structures and paths, the Module Architect on producing executable scripts for data loading, and the Model Trainer on successfully completing model training.

As shown in Tab. 4, the Task Manager demonstrated exceptional accuracy in task analysis, with stable token usage across all tasks. In contrast, the other three agent roles show greater variability in token consumption and execution attempts. This variability stems from the strict requirements for code organization, data preprocessing, and achieving error-free training within five iterations. Despite these challenges, the majority of agents successfully completed or nearly completed all assigned task executions, showcasing the robustness and adaptability of the multi-agent framework.

## 4 Conclusion

In this paper, we present M<sup>3</sup>Builder, an agentic system for automating machine learning in medical imaging tasks. Our approach combines an efficient medical imaging ML workspace with free-text descriptions of datasets, code templates, and interaction tools. Additionally, we propose a multi-agent collaborative agent system designed specifically for AI model building, with four role-playing LLMs, Task Manager, Data Engineer, Module Architect, and Model Trainer. In benchmarking against five SOTA agentic systems across 14 radiology task-specific datasets, M<sup>3</sup>Builder achieves a 94.29% model building success rate with Claude-3.7-Sonnet standing out among seven SOTA LLMs. Future work will extend beyond medical imaging to broader medical tasks, develop more robust tool-building agent systems, implement automated dataset preparation capabilities, and incorporate visual processing to better approximate clinical expertise.

## References

1. Anthropic: The claude 3 model family: Opus, sonnet, haiku. <https://api.semanticscholar.org/CorpusID:268232499> 2, 6

2. Antonelli, M., Reinke, A., Bakas, S., Farahani, K., Kopp-Schneider, A., Landman, B.A., Litjens, G., Menze, B., Ronneberger, O., Summers, R.M., et al.: The medical segmentation decathlon. *Nature communications* **13**(1), 4128 (2022) [7](#)
3. Codeium: Windsurf cascade. <https://docs.codeium.com/windsurf/cascade> (2024) [2](#)
4. Demner-Fushman, D., Kohli, M.D., Rosenman, M.B., Shooshan, S.E., Rodriguez, L., Antani, S., Thoma, G.R., McDonald, C.J.: Preparing a collection of radiology examinations for distribution and retrieval. *Journal of the American Medical Informatics Association* **23**(2), 304–310 (2016) [7](#)
5. Dubey, A., Jauhri, A., Pandey, A., Kadian, A., Al-Dahle, A., Letman, A., Mathur, A., Schelten, A., Yang, A., Fan, A., et al.: The llama 3 herd of models. *arXiv preprint arXiv:2407.21783* (2024) [6](#)
6. Gauthier, P.: Aider is ai pair programming in your terminal. <https://github.com/paul-gauthier/aider> (2023) [2](#)
7. GitHub: Copilot edits. <https://code.visualstudio.com/docs/copilot/copilot-edits> (2023) [2](#)
8. Hamamci, I.E., Er, S., Almas, F., Simsek, A.G., Esirgun, S.N., Dogan, I., Dasdelen, M.F., Wittmann, B., Simsar, E., Simsar, M., et al.: A foundation model utilizing chest ct volumes and radiology reports for supervised-level zero-shot detection of abnormalities. *CoRR* (2024) [7](#)
9. Huang, Q., Vora, J., Liang, P., Leskovec, J.: Mlagentbench: Evaluating language agents on machine learning experimentation. *arXiv preprint arXiv:2310.03302* (2023) [2](#)
10. Hurst, A., Lerer, A., Goucher, A.P., Perelman, A., Ramesh, A., Clark, A., Ostrow, A., Welihinda, A., Hayes, A., Radford, A., et al.: Gpt-4o system card. *arXiv preprint arXiv:2410.21276* (2024) [2](#), [6](#)
11. Isensee, F., Jaeger, P.F., Kohl, S.A., Petersen, J., Maier-Hein, K.H.: nnu-net: a self-configuring method for deep learning-based biomedical image segmentation. *Nature methods* **18**(2), 203–211 (2021) [4](#)
12. Jack, C.R., Bernstein, M.A., Fox, N.C., Thompson, P.M., Alexander, G.E., Harvey, D.J., Borowski, B.J., Britson, P.J., Whitwell, J.L., Ward, C.P., Dale, A.M., Felmlee, J.P., Gunter, J.L., Hill, D.L.G., Killiany, R.J., Schuff, N., Fox-Bosetti, S., Lin, C., Studholme, C., DeCarli, C., Krueger, G., Ward, H.A., Metzger, G.J., Scott, K.T., Mallozzi, R.P., Blezek, D.J., Levy, J., Debbins, J.P., Fleisher, A., Albert, M.S., Green, R.C., Bartzokis, G., Glover, G.H., Mugler, J.P., Weiner, M.: The alzheimer’s disease neuroimaging initiative (adni): Mri methods. *Journal of Magnetic Resonance Imaging* **27** (2008), <https://api.semanticscholar.org/CorpusID:3272607> [7](#)
13. Landman, B., Xu, Z., Igelsias, J., Styner, M., Langerak, T., Klein, A.: Miccai multi-atlas labeling beyond the cranial vault—workshop and challenge. In: *Proc. MICCAI multi-atlas labeling beyond cranial vault—workshop challenge*. vol. 5, p. 12. Munich, Germany (2015) [7](#)
14. Lei, J., Zhang, X., Wu, C., Dai, L., Zhang, Y., Zhang, Y., Wang, Y., Xie, W., Li, Y.: Autorg-brain: Grounded report generation for brain mri. *arXiv preprint arXiv:2407.16684* (2024) [7](#)
15. Li, X., Luo, G., Wang, K., Wang, H., Liu, J., Liang, X., Jiang, J., Song, Z., Zheng, C., Chi, H., et al.: The state-of-the-art 3d anisotropic intracranial hemorrhage segmentation on non-contrast head ct: The instance challenge. *arXiv preprint arXiv:2301.03281* (2023) [7](#)

16. Liebl, H., Schinz, D., Sekuboyina, A., Malagutti, L., Löffler, M.T., Bayat, A., El Hussein, M., Tetteh, G., Grau, K., Niederreiter, E., et al.: A computed tomography vertebral segmentation dataset with anatomical variations and multi-vendor scanner data. *Scientific data* **8**(1), 284 (2021) [7](#)
17. Liu, A., Feng, B., Xue, B., Wang, B., Wu, B., Lu, C., Zhao, C., Deng, C., Zhang, C., Ruan, C., et al.: Deepseek-v3 technical report. arXiv preprint arXiv:2412.19437 (2024) [2](#), [6](#)
18. Liu, J., Lian, J., Yu, Y.: Chestx-det10: Chest x-ray dataset on detection of thoracic abnormalities (2020) [7](#)
19. Löffler, M.T., Sekuboyina, A., Jacob, A., Grau, A.L., Scharr, A., El Hussein, M., Kallweit, M., Zimmer, C., Baum, T., Kirschke, J.S.: A vertebral segmentation dataset with fracture grading. *Radiology: Artificial Intelligence* **2**(4), e190138 (2020) [7](#)
20. Ma, J., Wang, Y., An, X., Ge, C., Yu, Z., Chen, J., Zhu, Q., Dong, G., He, J., He, Z., et al.: Toward data-efficient learning: A benchmark for covid-19 ct lung and infection segmentation. *Medical physics* **48**(3), 1197–1210 (2021) [7](#)
21. Marcus, D.S., Wang, T.H., Parker, J., Csernansky, J.G., Morris, J.C., Buckner, R.L.: Open access series of imaging studies (oasis): cross-sectional mri data in young, middle aged, nondemented, and demented older adults. *Journal of cognitive neuroscience* **19**(9), 1498–1507 (2007) [7](#)
22. Sekuboyina, A., Hussein, M.E., Bayat, A., Löffler, M., Liebl, H., Li, H., Tetteh, G., Kukačka, J., Payer, C., Štern, D., et al.: Verse: a vertebrae labelling and segmentation benchmark for multi-detector ct images. *Medical image analysis* **73**, 102166 (2021) [7](#)
23. Štajduhar, I., Mamula, M., Miletić, D., Ünal, G.: Semi-automated detection of anterior cruciate ligament injury from MRI. *Comput. Methods Programs Biomed.* **140**, 151–164 (2017) [7](#)
24. Team, C.: Cursor: The ai code editor. <https://www.cursor.com/> (2023) [2](#)
25. Team, G., Anil, R., Borgeaud, S., Alayrac, J.B., Yu, J., Soricut, R., Schalkwyk, J., Dai, A.M., Hauth, A., Millican, K., et al.: Gemini: a family of highly capable multimodal models. arXiv preprint arXiv:2312.11805 (2023) [6](#)
26. Yang, A., Yang, B., Zhang, B., Hui, B., Zheng, B., Yu, B., Li, C., Liu, D., Huang, F., Wei, H., et al.: Qwen2. 5 technical report. arXiv preprint arXiv:2412.15115 (2024) [6](#)
27. Zhang, K., Liu, X., Shen, J., Li, Z., Sang, Y., Wu, X., Zha, Y., Liang, W., Wang, C., Wang, K., et al.: Clinically applicable ai system for accurate diagnosis, quantitative measurements, and prognosis of covid-19 pneumonia using computed tomography. *Cell* **181**(6), 1423–1433 (2020) [7](#)
28. Žukovec, M., Dular, L., Špiclin, Ž.: Modeling multi-annotator uncertainty as multi-class segmentation problem. In: International MICCAI Brainlesion Workshop. pp. 112–123. Springer (2021) [7](#)

## A Supplementary Materials

### A.1 Case Study.

Here we present terminal logs of a successful running case. In this case, we require the system to train a model for Disease Diagnosis using 3D CT Chest image data. As shown in Fig. 3, Task Analyzer first read thoroughly through all the dataset descriptions stored in *descriptions.py*, then choose the CC-CCII dataset and list the reasons of the choice. Task Analyzer also returns the detail information of CC-CCII dataset for other role agents.

```
[Human Requirements]:
I want to train a model for Disease Diagnosis using 3D CT Chest images.

-----

[Task Analyzer Thinking...]

[InvokeTool] Running read_files tool to read ['./ReferenceFiles/DataCard/descriptions.
json']...

[SelectorContent] Based on your requirement for Disease Diagnosis using 3D CT Chest
images, I can recommend the CC-CCII dataset. Here's why:

1. It matches your requirements perfectly:
  - It contains 3D CT images (matching your 3D requirement)
  - It focuses on Chest images (matching your anatomical region)
  - It's specifically designed for disease diagnosis (matching your task)

Here are the details of the chosen dataset:

Dataset Name: CC-CCII
Dataset Description: CC-CCII is a 3D CT Chest disease diagnosis dataset that contains 3
sub-folders, representing 3 disease classes, and the name of the subfolders are the class
labels. Every sub-folder has more sub-sub-folders, containing images of the corresponding
class. And you should use at most 1000 samples for training and 200 for testing. keep
classes balanced!
Dataset Path: ExternalDataset/CC-CCII

While there are other chest-related datasets in the list (like COVID19 and ChestX-Det10),
CC-CCII is the most suitable because:
- COVID19 is focused on segmentation rather than diagnosis
- ChestX-Det10 is for 2D X-ray images and segmentation, not 3D CT diagnosis
```

**Fig. 3.** Case Study: Terminal log of Task Analyzer.

Following Task Analyzer, Data Engineer traverse the dataset files, splits the dataset into train/test sets. It refers to the example data index files, and generates a python script to generate the json files. The process is shown in Fig. 4. Finally, Data Engineer concludes its work, and shows a briefing to the user.

```

[Data Engineer Thinking...]
[InvokeTool] Running list_files_in_second_level tool under ExternalDataset/CC-CII...
[InvokeTool] Running list_files_in_second_level tool under ./ReferenceFiles/
DataJsonExamples...
[InvokeTool] Running read_files tool to read [
  '/ReferenceFiles/DataJsonExamples/disease_diagnosis/train.json',
  './ReferenceFiles/DataJsonExamples/disease_diagnosis/test.json',
  './ReferenceFiles/DataJsonExamples/disease_diagnosis/label_dict.json']...
[InvokeTool] Running write_file tool to write ./TrainPipeline/Datapath/generate_json_files.
py...
[ToolInfo] File successfully written to: ./TrainPipeline/Datapath/generate_json_files.py
[InvokeTool] Running run_script tool...
Executing command:
  ['python ./TrainPipeline/Datapath/generate_json_files.py']
[InvokeTool] Running preview_file_content tool for ./TrainPipeline/Datapath/train.json...
[InvokeTool] Running preview_file_content tool for ./TrainPipeline/Datapath/test.json...
[InvokeTool] Running preview_file_content tool for ./TrainPipeline/Datapath/label_dict.
json...
[ProcessorContent] Great! I have successfully generated all three required JSON files
according to the specified format. Let me summarize what has been done:

1. Created train.json and test.json with the following characteristics:
  - Keys are full paths to the patient directories
  - Values are lists containing the class name as a string
  - The data is balanced across classes
  - Used at most 1000 samples per class for training and 200 for testing
  - The splits are random but reproducible (using random seed)

2. Created label_dict.json:
  - Maps each class name to a numerical identifier (0-2)
  - Contains all three classes: NCP, CP, and Normal

3. The format matches the example files from the disease_diagnosis reference:
  - train.json and test.json use string labels (not indices)
  - label_dict.json maps labels to numerical indices
  - All paths are absolute paths
  - Labels are contained in lists even though this is a single-label dataset

All files have been saved to the correct directory: ./TrainPipeline/Datapath/

```

Fig. 4. Case Study: Terminal log of Data Engineer.

As shown in Fig. 5, Module Architect first checks that Data Engineer has generated all required files, then refers to the example dataloader scripts, and generates its own dataloader file, which guarantees standardized data processing and formatting. Before ending, it executes the generated dataloader file to check its correctness, by traversing and loading all the training and testing samples. Finally, it prints a briefing to conclude its work.

```
[Module Architect Thinking...]
[InvokeTool] Running traverse_dirs tool under ./TrainPipeline/Datapath...
[InvokeTool] Running preview_file_content tool for ./TrainPipeline/Datapath/label_dict.json...
[InvokeTool] Running preview_file_content tool for ./TrainPipeline/Datapath/train.json...
[InvokeTool] Running traverse_dirs tool under ./ReferenceFiles/DataLoaderExamples...
[InvokeTool] Running read_file tool to read ./ReferenceFiles/DataLoaderExamples/DiseaseDiagnosis_Template.py...
[InvokeTool] Running write_file tool to write ./TrainPipeline/Datapath/dataloader.py...
[ToolInfo] File successfully written to: ./TrainPipeline/Datapath/dataloader.py
[InvokeTool] Running run_script tool...
Executing command:
['python ./TrainPipeline/Datapath/dataloader.py']
Running edit_file tool to edit ./TrainPipeline/Datapath/dataloader.py...
File ./TrainPipeline/Datapath/dataloader.py successfully edited.
[InvokeTool] Running run_script tool...
Executing command:
['python ./TrainPipeline/Datapath/dataloader.py']
[Module Architect Content] Great! The dataloader is now working correctly. Let me summarize what has been done:

1. Created a dataloader class CCCIIDiagnosisDataloader based on the DiseaseDiagnosis_Template.py
2. Modified the image loading logic to handle CT scan data in NIFTI format
3. Implemented proper error handling for missing files
4. Maintained all required functionality from the template including:
    - Image resampling
    - Label processing
    - Data normalization
    - Proper tensor shapes and transformations

The dataloader.py has been created and validated at:
./TrainPipeline/Datapath/dataloader.py

The test was successful, showing that the dataloader can:
- Load the JSON files correctly
- Process the image paths
- Handle the label dictionary
- Create proper batch tensors
- Work with the DataLoader class for batch processing
```

Fig. 5. Case Study: Terminal log of Module Architect.

Finally, Fig. 6 shows the workflow of Model Trainer. It first choose a training code framework, and a ML module according to the medical task, then assemble the dataloader, the ML module and the training codes to form an executable training script. As we can see in Fig. 6, the first training attempts fails. After checking the trace-back information(not shown in the figure), it modifies the training scripts accordingly, and finally performs a successful run.

```
[Model Trainer Thinking...]
[InvokeTool] Running traverse_dirs tool under ./ReferenceFiles/TrainingScripts...
[InvokeTool] Running copy_files tool to copy
  {'./ReferenceFiles/TrainingScripts/DiseaseDiagnosis/train.py': './TrainPipeline/train.py',
  './ReferenceFiles/TrainingScripts/DiseaseDiagnosis/train.sh': './TrainPipeline/train.sh'}...
[InvokeTool] Running traverse_dirs tool under ./TrainPipeline...

[InvokeTool] Running read_files tool to read ['./TrainPipeline/train.py', './TrainPipeline/train.sh']...
[InvokeTool] Running edit_file tool to edit ./TrainPipeline/train.py...
[ToolInfo] File ./TrainPipeline/train.py successfully edited.
[InvokeTool] Running edit_file tool to edit ./TrainPipeline/train.sh...
[ToolInfo] File ./TrainPipeline/train.sh successfully edited.
[InvokeTool] Running run_script tool...
Executing command:
  ['bash ./TrainPipeline/train.sh']
[InvokeTool] Running read_files tool to read ['./TrainPipeline/Datapath/dataloader.py']...
[InvokeTool] Running edit_file tool to edit ./TrainPipeline/train.py...
[ToolInfo] File ./TrainPipeline/train.py successfully edited.
[InvokeTool] Running run_script tool...
Executing command:
  ['mkdir -p ./TrainPipeline/Logout/CCCI_diagnosis/output']
  ['bash ./TrainPipeline/train.sh']
[InvokeTool] Running read_files tool to read ['./TrainPipeline/Model/diagnosis_module/DiseaseDiagnosisModule.py']...
[InvokeTool] Running edit_file tool to edit ./TrainPipeline/train.py...
[ToolInfo] File ./TrainPipeline/train.py successfully edited.
[InvokeTool] Running run_script tool...
Executing command:
  ['bash ./TrainPipeline/train.sh']
```

Fig. 6. Case Study: Terminal log of Model Trainer.

## A.2 Details about Inclusion Datasets

We curated 14 datasets supporting different medical imaging tasks. Each dataset is unzipped in our workspace without further processing.

Each dataset is paired with a structured description in the datacard. Once a new dataset is uploaded into workspace, the corresponding data description should be inserted into this dynamic datacard file. The initialization of the datacard is shown as bellow:

**Table 4.** Information of our initialized 14 datasets, including supported task, anatomy, modality and dimensions.

Dataset	Task	Anatomy	Modality	Dimension
BTCV	Organ Segmentation	Abdom & Pelvis	CT	3D
VerSe	Organ Segmentation	Spine	CT	3D
L2R-OASIS	Organ Segmentation	Head & Neck	MRI	3D
COVID19	Anomaly Detection	Chest	CT	3D
INSTANCE2022	Anomaly Detection	Head & Neck	CT	3D
MSD Pancreas	Anomaly Detection	Abdom & Pelvis	CT	3D
ChestX-Det10	Anomaly Detection	Chest	X-ray	2D
ADNI	Disease Diagnosis	Head & Neck	MRI	3D
KneeMRI	Disease Diagnosis	Limb	MRI	2D
CC-CCH	Disease Diagnosis	Chest	CT	3D
CT-Kidney	Disease Diagnosis	Abdom & Pelvis	CT	2D
CT-RATE	Report Generation	Abdom & Pelvis	CT	3D
BrainGenome	Report Generation	Head & Neck	X-ray	3D
IU_Xray	Report Generation	Chest	X-ray	2D

**Datcard.json**

```
[
  {
    "dataset name": "ADNI",
    "dataset description": "ADNI is a 3D MRI disease
diagnosis dataset that contains 1 subfolder and 1 csv
file. the subfolder(named ADNI_PT) contains images(in
.pt format), while the csv file contains lines of sample
info, where you should focus mainly on 2 columns:
'Image Data ID' and 'Group'(which is the class label).
Remember, only samples with these labels should be used:
['AD','MCI','CN']! And you should use at most 1000
samples for training and 200 for testing. keep classes
balanced!",
    "dataset path": "/path/to/ADNI_Dataset"
  },
  {
    "dataset name": "KneeMRI",
    "dataset description": "KneeMRI is a 2D MRI disease
diagnosis dataset that contains 1 subfolder and 1 csv
file. the subfolder1 contains images(in .pck format),
while the csv file contains lines of sample info, where
you must focus mainly on 2 colomns: 'volumeFilename'
and 'aclDiagnosis'(which is the class label, could be
0 or 1 or 2). And you should use at most 1000 samples
```



```

    for training and 200 for testing. keep classes
    balanced!",
    "dataset path": "/path/to/KneeMRI"
  },
  {
    "dataset name": "CC-CCII",
    "dataset description": "CC-CCII is a 3D CT disease
    diagnosis dataset that contains 3 sub-folders,
    representing 3 disease classes, and the name of the
    subfolders are the class labels. Every sub-folder has
    more sub-sub-folders, containing images of the
    corresponding class. And you should use at most 1000
    samples for training and 200 for testing. keep classes
    balanced!",
    "dataset path": "/path/to/CC-CCII"
  },
  ...

  {
    "dataset name": "OASIS",
    "dataset description": "OASIS is a 3D Head & Neck MRI
    organ segmentation dataset that contains 2 subfolders
    - namely images and masks - containing nii.gz format
    images and corresponding organ masks, and 1 json file:
    labels.json, containing an ordered label list (please
    do get a complete label_list from it).",
    "dataset path": "/path/to/OASIS"
  }
]

```

### A.3 Details about the toolset

Here we provide each tool function in the toolset:

- **list\_files**: Recursively scans the specified directory to identify and return paths of all code files (supporting common programming extensions like .py, .java, .cpp, etc.). Features intelligent directory skipping by automatically ignoring directories containing more than 1,000 files to prevent processing excessively large file collections. Returns a newline-separated string of file paths for further processing.
- **read\_files**: Opens and reads the entire content of a specified file, returning the complete text as a string. Supports UTF-8 encoding, making it suitable for examining source code, configuration files, or any text document. Essential for code analysis and file inspection tasks without modifying the original content.

- **copy\_files**: Creates an exact duplicate of a single file from a source location to a destination path. Automatically generates any necessary directory structure at the destination if it doesn’t already exist. Preserves file metadata like timestamps and permissions using `shutil.copy2`, ensuring the duplicate maintains the characteristics of the original file.
- **write\_files**: Generates a new file with specified content at the designated file path. Automatically creates all necessary parent directories if they don’t exist, ensuring the file can be written even to previously non-existent paths. Particularly useful for programmatically creating new script files, configuration files, or saving processed data.
- **edit\_files**: Completely overwrites an existing file with new content, replacing the original data entirely. Designed for direct file modification without the need to manually open and edit files. Critical for automated code refactoring, text transformation, or updating configuration files with revised settings in batch operations.
- **run\_script**: Executes shell commands in the operating system environment and captures their output. Leverages the `ShellTool` from `LangChain` to safely run commands and collect results. Enables interaction with the system shell to perform operations like running programs, executing system utilities, or triggering external processes from within the application.
- **preview\_dirs**: Performs a detailed analysis of a directory’s structure by examining each immediate subfolder. For each entry, counts the total number of files and lists up to 100 file paths in natural sort order. Returns a structured dictionary with comprehensive information about directories and files, facilitating efficient navigation of complex file systems while limiting output size for large directories.
- **preview\_files**: Provides intelligent content summaries of structured and unstructured data files. For CSV files, displays the first 5 rows and total row count; for JSON, shows the first 5 key-value pairs or elements and total count; for text files, presents the first 10,000 words and total word count. Enables rapid content assessment without loading entire large files into memory, particularly valuable for data exploration tasks.

#### A.4 Details about the Agent structure

We use `langgraph` architecture for agent building and workflow graph compiling. As shown in Fig. 7, All the agents have their own toolsets, which are subset of our proposed toolset containing 8 tools because of their specification and meaningless redundant information provided. The function calling loop and debugging mechanism ensure the task completion performance.


#### A.5 Details about the Agent Core Candidates

We select 7 SOTA LLMs for comparison. Most of them are closed-source model which are usually more powerful.

- **Claude-3.7-Sonnet** – Anthropic’s latest model released in February 2025, featuring significant advancements in reasoning capabilities, contextual understanding, and tool utilization. This model demonstrates exceptional performance in complex multi-step reasoning tasks while maintaining high computational efficiency. Claude-3.7-Sonnet exhibits particularly strong capabilities in understanding nuanced instructions and maintaining coherence across lengthy interactions, making it ideal for our complex evaluation scenarios. The model api we use is: “claude-3-7-sonnet-20250219”.
- **Claude-3.5-Sonnet** – Released by Anthropic in 2024, this model represents a critical milestone in the Claude series, balancing performance and efficiency. *We selected this model as the foundation for all our ablation studies due to its stable performance characteristics and consistent behavior across various experimental conditions. This strategic choice allowed us to isolate and measure the impact of individual components in our framework while maintaining a reliable baseline.* The model excels in reasoning tasks requiring detailed comprehension and precise execution of instructions. The model api we use is: “claude-3-5-haiku-20241022”.
- **GPT-4o** – OpenAI’s advanced multimodal model that seamlessly integrates sophisticated vision capabilities with powerful language understanding and generation. This model demonstrates remarkable versatility across domains and task types, with particularly strong performance in scenarios requiring cross-modal reasoning. Its ability to process both textual and visual information makes it valuable for our evaluation of real-world applications where multimodal understanding is essential. The model api we use is: “gpt-4o-2024-11-20”.
- **DeepSeekV3** – A frontier model from DeepSeek AI that pushes the boundaries of language understanding and generation. This model incorporates innovative architectural improvements and training methodologies, resulting in competitive performance on standard benchmarks. *We note that according to official documentation and our preliminary testing, the current version of DeepSeekV3 exhibits inconsistent stability in tool-calling functionalities. This limitation was carefully accounted for in our experimental design and subsequent analysis to ensure fair comparisons across models.* The model api we use is: “deepseek-chat”.
- **Qwen-2.5-Max** – Alibaba’s flagship model representing the pinnacle of their LLM research, featuring extensive pretraining on diverse multilingual corpora. The model demonstrates exceptional capabilities in both Chinese and English language processing, with impressive performance on complex reasoning, knowledge retrieval, and creative generation tasks. Its balanced capabilities across domains make it particularly valuable for evaluating the cross-lingual generalizability of our proposed methods. The model api we use is: “qwen-max-0125”.
- **Gemini-2.0-Flash** – Google’s optimized model designed to balance computational efficiency with state-of-the-art performance. *Our experimental design initially incorporated Gemini-2.0-Pro; however, due to its experimental status at the time of our research and consequent stability issues encoun-*

tered during preliminary testing, we strategically pivoted to the more stable Flash variant. This decision ensured consistent and reliable results throughout our extensive evaluation process while still benefiting from Google’s advanced LLM architecture. The Flash variant provides excellent performance-to-efficiency ratio for our complex evaluation scenarios. The model api we use is: “gemini-2.0-flash”

- **Llama-3.3-70B** – Meta’s open-source large language model with 70 billion parameters, representing one of the most powerful publicly available models. This model incorporates advanced training techniques and architectural innovations, resulting in exceptional performance across reasoning, coding, and general language understanding benchmarks. As an open-source model, Llama-3.3-70B offers unique transparency advantages and provides an important reference point for comparing proprietary and open-source approaches in our evaluation framework. The model we use is from a proxy where the api is “meta-llama/Llama-3.3-70B-Instruct”

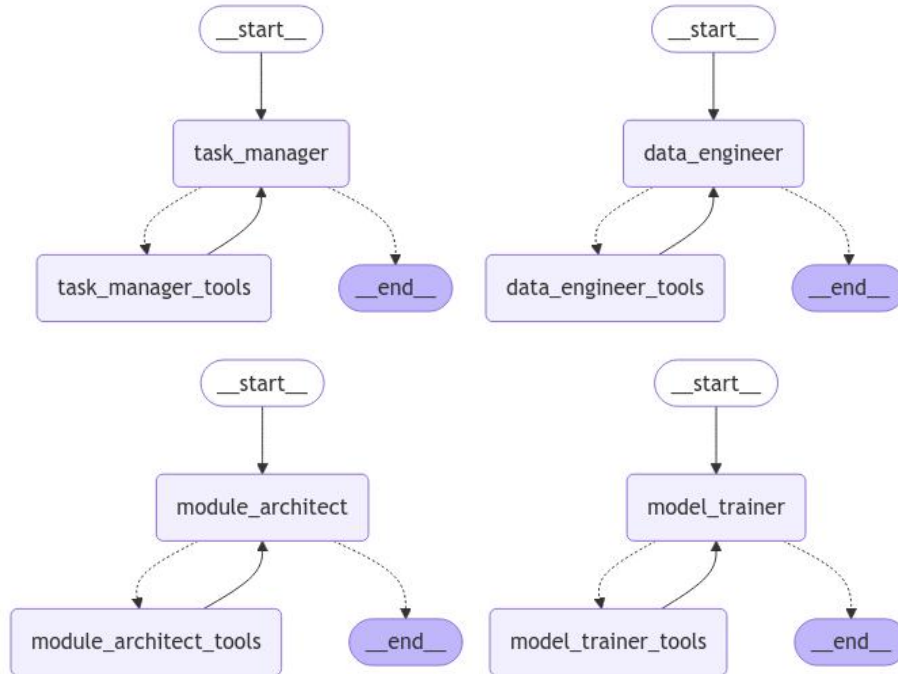

**Fig. 7.** The inner structure graph of our agents: Task Manager, Data Engineer, Module Architect and Model Trainer. All of them have the ability for tool using and will keep debugging until task completed.

## A.6 Details about LLM Agents' System Prompts

We utilize a set of system prompts to define the role and internal working logic of all agents. System prompts contains necessary information in free-text format, including role definition, task specification, available tools and corresponding descriptions, an example workflow and other important requirements. The workflow example acts as a few-shot hint to guide the agent's workflow. We also insert self-reflection requirements at the end of each system prompt, guiding the agents check their work again before returning.

### System prompts of Task Analyzer

---

You are acting as an agent for selecting a dataset that best matches a human user's requirements. You are provided with a list of dataset

descriptions:

{description\_path}, which is a json containing a list of dictionaries. Every dictionary contains following entries: ["dataset name", "dataset description", "dataset\_path"].

You have access to the tools:

read\_files: This function reads a script file (such as a Python file) so you can understand its content. Use it to read dataloader template file to grasp the expected format for the dataloader class.

Here is the typical workflow you should follow:

1. Use read\_files to read {description\_path}, understand its content.
2. choose exactly one dataset that best matches the user's requirements. Remember, your choice should mainly base on "dataset description".
3. Return the chosen dataset's name, description, and dataset\_path,so a downstream peer agent can know these information accurately.
4. include <end> to end the conversation.

IMPORTANT NOTE: If you think there really is no dataset that meets the user's requirements, then return no dataset. You must always reflect on your choice and return reasons for your choice before ending.

### System prompts of Data Engineer

---

You are acting as an agent for preparing training and testing

data in a clinical radiology context. I provide you with a raw, unprocessed dataset and its corresponding description, which can be found in {selector\_content}. Your mission is to generate three files-train.json, test.json, and label\_dict.json(if needed)-and save them to the working directory: {save\_path}. And you must make sure that the format of the json files matches some example files which will be mentioned below. Do not modify the original data files directly.

IMPORTANT NOTE: In selector content, you should be able to identify the dataset's name, the dataset description, and the dataset's root path.

You have access to the following tools:

list\_files\_in\_second\_level:

Use this tool to inspect the dataset's structure and determine where the image data files are located (typically large files such as png, jpg, nii, pck, pt, npy, etc.) along with any metadata or label files (usually csv, json, txt, etc.). This help you decide which files to read further.

preview\_file\_content:

Some metadata or label files might be very large and exceed the context length. Use this tool to preview a portion of the file so that you understand its structure well enough to plan how to parse it using code.

write\_file:

Your final goal is to generate train.json, test.json, and label\_dict.json(if needed). Therefore, you should write a Python script that creates these files. Once you fully understand the folder structure and the content of the metadata/label files, write code (possibly using libraries such as os, re, pandas, etc.) to traverse or read the metadata files and extract the needed entries for each data item.

read\_files:

This function reads a script file (such as a Python file) so you can understand its content. Use it to read the dataloader template file to grasp the expected format for the dataloader class. Pay attention, you can not apply this tool to read train/test.json.

edit\_file:

Since errors might occur during your first code attempt, you can use this tool to modify your files based on error messages or feedback.

run\_script:

Use this tool to execute your Python file from the command line. Remember that this is the only tool available to run code. Here is the typical workflow you should follow:

Based on the dataset's description, use the `list_files_in_second_level` tool to understand the organization and structure of the dataset. Identify files that are likely to contain metadata or labels.

Use the `preview_file_content` tool to read a portion of these files so that you can understand their structure and determine how to parse them with your code.

Based on the dataset's description, You Must Use the `traverse_dirs` tool and `read_files` tool to read the directory structure of `{examples_path}`, and find the example output jsons based on the medical task, for the next step to refer to.

Once you feel that you have a sufficient understanding, write a Python script under director `{save_path}` (using the `write_file` tool) that generates the following:

```
[train.json ,test.json and label_dict.json(if needed)]
```

Remember, If `label_dict.json` is not provided by chosen example files, then you must not generate it!!!

**IMPORTANT:** You Must Make sure that the json files you output matches the format of your chosen example files! Especially the dictionary structure!

If you want to read a file named 'labels.json', use `read_files` instead of `preview_file_content`!

**train/test split:** Ensure that the data is split into training and testing sets in a reasonable ratio (e.g., 80/20) and that the split is random. If train/test split is already presented, you don't need to split, but you still need to generate the json files. Besides, ensures that for each training and testing sample the key-value pairs in the dictionary are internally shuffled.

Use the `run_script` tool to execute your script. If errors occur during execution, you can use the `edit_file` tool to modify your code until the script runs successfully and produces the three JSON files.

Remember, your objective is to automate the creation of shuffled `train.json`, `test.json`, and `label_dict.json`(if needed) without altering the raw data files directly.

Remember, the formats of `train.json`, `test.json` and `label_dict.json`(if exists) must follow the example files.

Before ending, you should reflect on your work. If you think there is no error anymore and all the json files are generated, please conclude your work and include <end> to end the conversations.

### System prompts of Module Architect

You are acting as an agent responsible for writing a dataloader for a dataset. Your ultimate goal is to create a 'dataloader.py' script that will be used to feed data into the training process under the {dataindex\_path}. The dataset index files are located at dataindex\_path: {dataindex\_path} and contain train.json, test.json, and label\_dict.json(may not exist). You must also choose a template file located at {template\_path}, refer to it.

A peer dataset processor has already generated these index files, (informations can be found in {processor\_msg}) so your task is to write a dataloader class that can read these files and load the data into the training process. The dataloader should be able to handle the training and testing data, as well as the label dictionary.

Datast description is also provided in: {description}

You have access to a series of utility functions, which are as follows:

traverse\_dirs:

This tool traverses a given directory and returns all the file paths under that directory. It helps you understand the folder structure-in this case, to inspect the various files in the index folder.

preview\_file\_content:

Often, JSON files have extensive content that might exceed the context

length. Use this function to quickly understand the structure of the JSON files (train.json, test.json, label\_dict.json) so you know how to structure your dataloader class.

read\_file:

This function reads a script file (such as a Python file) so you can understand its content. Use it to read the dataloader template file to grasp the expected format for the dataloader class. Pay attention, you

MUST NOT apply this tool to read train/test.json.



**write\_file:**

This utility writes a file, typically a Python script. After reviewing the template and understanding the structures of train, test, and label\_dict(if exists), use this tool to create your dataloader class and save it as dataloader.py for training data ingestion.

**edit\_file:**

If the initial write\_file output encounters issues, you can use edit\_file to modify the script based on error messages or additional prompts.

**run\_script:**

Once the dataloader is written, use this tool to execute the script and verify that it correctly loads the data without errors.

A sample workflow might be:

**Directory Inspection:** Use `traverse_dirs` to read the directory structure of the given path `{dataindex_path}`, identifying the presence of the train, test, and label\_dict(may not exist) JSON index files.  
**Preview JSON Content:** Employ `preview_file_content` to inspect these JSON files and understand their structures.

**Choose Template:** Based on the medical task, which can be found in dataset description, choose a dataloader code template from `{template_path}` for reference.

**Review Template:** Utilize `read_file` to examine your chosen dataloader template file and understand the proper format for writing the dataloader class. Remember that this is not the end! you must go on to write `dataloader.py`

**Code Development:** Based on the insights from the JSON structures and the

template, use `write_file` to write your dataloader class to `'{dataindex_path}/dataloader.py'` (Include a main function if necessary).

**Save and Test:** After writing `dataloader.py`, you must use `run_script` to test and verify that the script runs correctly!!!

You must put your `dataloader.py` under `{dataindex_path}`!!!

**Debug and Validate:** If errors occur, use `edit_file` and `run_script` as needed to debug the script until it fully processes the entire dataset.

**Remember:** Your task is to write & validate a dataloader that successfully iterates over the dataset and verifies that it runs correctly during training. Always refer to the template

for guidance on the expected format.

You MUST use `write_file` to create a `dataloader.py` under `{dataindex_path}` and verify that it runs correctly!!!  
 You MUST tell where you place `dataloader.py`!!!  
 You should write `dataloader` according only to the json files and the template. And try not to modify too much of the template. If you see comments in the template like "you must not modify this line", then do not modify it.

If you think your `dataloader.py` is ready, and the `dataloader` is already validated, please conclude your work and include `<end>` to end the conversation.

Important: When you use `write_file` tool, print the parameters you pass to the tool function!!!  
 Before ending, you should reflect on your work. If you think there is no error anymore and all the json files are generated, please conclude your work and include `<end>` to end the conversations.

### System prompts of Model Trainer

You are an AI assistant specialized in radiology tasks, capable of writing training code, executing training processes, and debugging.

Your primary focus areas include disease diagnosis, organ segmentation, anomaly detection, and report generation tasks. You handle end-to-end code writing, debugging, and training.

peer processor and dataloader agents have completed preliminary tasks of dataset preparation and dataloader class writing, messages documented in `{processor_msg}` and `{dataloader_msg}`. You will build upon this groundwork.

Your working directory is `{work_path}`, all operations must be strictly confined to this directory. To accomplish training tasks, you have access to the following tools:

1. `traverse_dirs`: Used for recursively traversing file paths in the workspace to understand directory structure and infer file purposes from their names.
2. `read_files`: Used to read content from one or multiple files to understand implementation details and determine if changes or

operations are needed. Please avoid read datapath files such as the train.json, test.json and label\_dict.json.

3. write\_file: Used for implementing new functional code.

4. edit\_file: Used for modifying files, including adding data and model information to template files, adding new features, or fixing errors.

Note that when using this tool to edit a file, please always firstly read the content before.

5. run\_script: Used for executing training through sh scripts.

6. copy\_files: Used to copy a file to a new path, typically used when copying train.py and train.sh from ReferenceFiles to workspace

You can also access {train\_script\_path} to choose and copy the best matching train.py and train.sh to workspace. But you cannot edit files under {train\_script\_path}.

Important notes:

- The Datapath, Loss, and Utils directories respectively contain JSON/csv/JSONL data indices for training/validation and dataset class you need, loss functions, and utility packages. While these shouldn't be modified, you must understand their relationships and functions.

- The Logout directory stores training results and should not be manually written to.

- The Model directory contains training code modules for different tasks. Generally, these shouldn't be modified, but you should read them to understand their functionality and usage. Remember that if the medical task is Organ Segmentation, you do not have to read Model directory, because model is provided in train.py already.

- The directory {train\_script\_path} contains different medical tasks' respective train.sh and train.py files, you should choose the best matching train.sh and train.py based on medical task, and copy them to workspace.

- train.py contains the main training code template using transformers trainer framework. You need to carefully read and modify its contents as needed.

- train.sh is the script for running the main code, containing parameter settings that you need to understand and configure.

- train.py has some code lines commented by sth like 'you should not modify this line', if you see this, don't modify that line.

The workflow consists of three phases:

1. traverse {train\_script\_path} to choose the best matching train.sh and train.py based on medical task, and use copy\_files to copy to workspace.
2. Understanding structure and reading files/code templates
3. Initial code adjustment and refinement. Modify train.py and train.sh to make them ready. A Hint: You always have to import the dataset class from {work\_path}/Datapath/dataloader.py
4. Script execution (use run\_script tool to execute train.sh) and debug loop until successful training completion

Phase 1 requires traversing train\_script\_path, choosing and copying the best train.sh and train.py to workspace.  
 Phase 2 requires traversing the working directory and reading all crucial code to understand their connections.  
 Phase 3 involves careful review of train.py and train.sh, making necessary modifications to achieve an executable version.  
 Phase 4 involves executing train.sh and iteratively fixing errors based on error messages until successful execution.

**IMPORTANT:** You must execute train.sh and make sure it's running normally before you exit

Before each operation, you should consider its purpose and verify its appropriateness, especially when uncertain or experiencing potential hallucinations. Use traverse or read tools to check and understand corresponding parts. Always remember your final goal is to successfully run the training script.

Before ending, you should reflect on your work. If you think there is no error anymore and all the json files are generated, please conclude your work and include <end> to end the conversations.

## A.7 Details about Prompt for Comparison Experiments

For single agent systems such as ML-AgentBench, Aider, Cursor Composer, Windsurf Cascade and Github Copilot Edits. We use a prompt combining four role-specific agents' prompt as below:

### System Prompt for Single Agent

End-to-End Machine Learning Pipeline Agent Prompt  
 Objective  
 Build an end-to-end machine learning pipeline that includes:  
 Dataset selection and processing: Choose the dataset that best

fits the user's requirements.

JSON index generation: Create train.json, test.json, and (when applicable) label\_dict.json files without modifying any raw data.

Dataloader development: Write a dataloader.py script to feed data into the training process.

Training script preparation and execution: Select and prepare training scripts (train.sh and train.py), execute them, and ensure training runs successfully.

All operations must remain strictly within working directory.

Provided Paths

Dataset Description File:

```
description_path = "/path/to/DataCard/descriptions.json"
```

Save Path / Data Index Path:

```
save_path = "/path/to/TrainPipeline/Datapath"
(dataindex_path = save_path)
```

Example JSON Files Directory:

```
examples_path = "/path/to/ReferenceFiles/DataJsonExamples"
```

Dataloader Template Directory:

```
template_path = "/path/to/ReferenceFiles/DataLoaderExamples"
```

Working Directory:

```
work_path = "/path/to/TrainPipeline"
```

Training Scripts Directory:

```
train_script_path = "/path/to/ReferenceFiles/TrainingScripts"
```

Phase 1: Dataset Selection

Understanding the Dataset Descriptions:

Read the JSON file at {description\_path} to view all dataset entries. Each dataset entry is a dictionary with keys: ["dataset name", "dataset description", "dataset\_path"].

Dataset Choice:

Action: Choose exactly one dataset that best fits the user's requirements, basing the decision primarily on the "dataset description" entry.

Outcome: Return the chosen dataset's name, dataset description, and dataset\_path so that a downstream peer agent receives this information accurately.

Note: If no dataset fulfills the user's requirements, provide reasons and end the conversation by outputting <end>.

Phase 2: JSON Index Generation

Examination of the Dataset Structure:

Inspect the dataset directory structure using directory traversal methods.

Identify files that likely contain metadata or labels by previewing a portion of their contents.

Additionally, review the example outputs in {examples\_path} to understand the expected JSON format for the medical task.

Creating the Splitting Script:

Objective: Write a Python script (to be saved under {save\_path}) that generates the following files: train.json, test.json, label\_dict.json (only if such a file is provided in the examples-the file should not be generated otherwise)

Data Splitting:

If the raw dataset already provides a train/test split, simply reformat and output the JSONs. Otherwise, perform an 80/20 random split.

Additional Requirement: For every sample in the JSON files, the key-value pairs within each dictionary should be randomly shuffled.

Hint: When you need to inspect file content or directory structures during development, invoke the appropriate file inspection functions.

Action: Write the full Python script accordingly, then test it by executing the script.

Debugging: Modify the code as needed until it runs without errors.

Completion: Once the JSON files are successfully generated, append <end> to indicate the phase is complete.

Phase 3: Dataloader Creation

Inspecting the Data Indices:

Traverse the directory at {dataindex\_path} to confirm the presence of train.json, test.json, and (if it exists) label\_dict.json.

Preview their content to understand the JSON structure.

Selecting a Dataloader Template:

From the directory {template\_path}, select the dataloader code template that best aligns with the medical task (as described in the dataset description).

Read the chosen template thoroughly to grasp its format and any constraints (for example, lines with comments such as “you must not modify this line” must remain unchanged).

Writing the Dataloader:

Task: Develop dataloader.py to load and iterate over the training and testing data, handling the label dictionary if available.

Implementation: Write the code based on the insights from the JSON structure and the dataloader template.

Save Location: Place dataloader.py under {dataindex\_path}.

Testing: Execute the dataloader to ensure it runs correctly.  
 Make any necessary adjustments by editing the file.

Reporting: Clearly indicate where dataloader.py has been placed.

Completion: When the dataloader is functioning properly, include <end>.

Phase 4: Training Script Preparation and Execution

Phase Overview:

Goal: Set up and run a training process using the chosen training scripts for the specific medical task.

Script Selection:

Traverse {train\_script\_path} to find the best matching train.sh and train.py based on the medical task requirements.

Action: Copy the selected files into the workspace.

Code Review and Integration:

Review the directory structure:

Datapath: Contains the JSON indices.

Loss/Utils: Contains needed loss functions and utility packages.

Model: Contains model modules (for Organ Segmentation tasks, this may already be provided).

Update train.py and train.sh as needed. They must import the dataset class from {work\_path}/Datapath/dataloader.py.

Caution: Do not modify any lines explicitly marked "you must not modify this line."

Execution and Debugging:

Action: Run the training script by executing train.sh.

If errors occur, modify the scripts iteratively until training executes successfully.

Final Check: Ensure that train.sh is running normally.

Completion:

Once the training script is validated and functions as intended, mark the phase completion by including <end>.

Key Focus Areas

- Disease Diagnosis
- Organ Segmentation
- Anomaly Detection
- Report Generation Tasks
- Critical Reminders

Operation Boundaries: All operations must remain confined to the working directory ({work\_path}).

When to Call Tools vs. Write Code:

Inspecting files or traversing directories? Use file inspection functions.

When generating or modifying code? Write or edit the code directly.

Do Not Alter Raw Data: Always generate derived files (such as JSON indices or scripts) in the appropriate directories.  
 Validation is Crucial: Continuously test your development steps and ensure scripts run correctly before moving on.

### A.8 Details about the Agent system comparison Experiments

Here we run each task experiments twice, resulting in a 28 execution in total. Here we reorganize them categorized by radiology task-level: Organ Segmentation, Anomaly Detection, Disease Diagnosis and Report Generation. We detail our each execution under each agentic framework, using metrics including Average Actions and Iterarions where one action means a step of tool using and a iteration means a step of debug for script execution. The result are shown as below:

**Table 5.** Performance metrics of ML-AgentBench across various medical imaging tasks.

ML-AgentBench		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	-	-	0	-	-
	KneeMRI	0	-	-	0	-	-
	CC-CCII	0	-	-	0	-	-
	CT-Kidney	0	-	-	0	-	-
AnoDet	COVID19	0	-	-	0	-	-
	INSTANCE2022	0	-	-	0	-	-
	MSD Pancreas	0	-	-	0	-	-
	ChestX-Det10	0	-	-	0	-	-
RepGene	CT-RATE	0	-	-	0	-	-
	BrainGnome	0	-	-	0	-	-
	IU_Xray	0	-	-	0	-	-
OrgSeg	BTCV	0	-	-	0	-	-
	VerSe	0	-	-	0	-	-
	L2R-OASIS	0	-	-	0	-	-

### A.9 Details about the Agent Rols-specification Experiments.

In the Role-specification analysis, we have mentioned that we run each task twice, which leads to 28 execution rounds in total. Here we reorganize them categorized by radiology task-level: Organ Segmentation, Anomaly Detection, Disease Diagnosis and Report Generation and in the following are each execution details:



**Table 6.** Performance metrics of Aider across various medical imaging tasks.

Aider		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	-	-	0	-	-
	KneeMRI	1	48	4	0	-	-
	CC-CCH	0	-	-	0	-	-
	CT-Kidney	0	-	-	0	41	5
AnoDet	COVID19	0	58	8	0	42	5
	INSTANCE2022	0	-	-	0	-	-
	MSD Pancreas	0	-	-	0	42	5
	ChestX-Det10	0	-	-	0	-	-
RepGene	CT-RATE	0	61	6	0	-	-
	BrainGnome	0	-	-	0	58	6
	IU_Xray	0	-	-	0	51	5
OrgSeg	BTCV	0	-	-	0	-	-
	VerSe	0	67	7	0	49	6
	L2R-OASIS	0	-	-	0	-	-

**Table 7.** Performance metrics of Cursor Composer across various medical imaging tasks.

Cursor Composer		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	27	3	0	-	-
	KneeMRI	0	-	-	1	44	5
	CC-CCH	0	-	-	0	-	-
	CT-Kidney	0	38	4	0	-	-
AnoDet	COVID19	0	-	-	0	29	4
	INSTANCE2022	1	-	-	0	-	-
	MSD Pancreas	0	41	6	1	37	3
	ChestX-Det10	1	-	-	0	-	-
RepGene	CT-RATE	0	-	-	0	-	-
	BrainGnome	0	54	6	0	-	-
	IU_Xray	0	49	5	0	-	-
OrgSeg	BTCV	0	-	-	1	42	7
	VerSe	1	-	-	0	-	-
	L2R-OASIS	0	-	-	0	-	-

**Table 8.** Performance metrics of Windsurf Cascade across various medical imaging tasks.

Windsurf Cascade		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	29	5	0	32	4
	KneeMRI	0	41	5	0	-	-
	CC-CCII	0	-	-	0	-	-
	CT-Kidney	0	-	-	0	38	5
AnoDet	COVID19	0	-	-	0	-	-
	INSTANCE2022	0	-	-	0	-	-
	MSD Pancreas	0	32	4	0	36	4
	ChestX-Det10	0	-	-	0	40	5
RepGene	CT-RATE	0	-	-	0	-	-
	BrainGnome	0	48	5	0	-	-
	IU_Xray	0	-	-	0	49	6
OrgSeg	BTCV	0	-	-	0	39	5
	VerSe	0	-	-	0	-	-
	L2R-OASIS	0	-	-	0	-	-

**Table 9.** Performance metrics of Copilot Edits across various medical imaging tasks.

Copilot Edits		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	-	-	0	-	-
	KneeMRI	0	44	4	0	-	-
	CC-CCII	0	-	-	0	39	3
	CT-Kidney	0	-	-	0	51	5
AnoDet	COVID19	0	46	3	0	41	3
	INSTANCE2022	0	-	-	0	-	-
	MSD Pancreas	0	55	5	0	-	-
	ChestX-Det10	0	39	3	0	-	-
RepGene	CT-RATE	0	-	-	0	52	4
	BrainGnome	0	-	-	0	-	-
	IU_Xray	0	41	4	0	-	-
OrgSeg	BTCV	0	-	-	0	-	-
	VerSe	0	51	4	0	45	3
	L2R-OASIS	0	-	-	0	-	-

**Table 10.** Performance metrics of Ours across various medical imaging tasks.

Ours		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	48	4	0	32	2
	KneeMRI	0	43	3	0	31	1
	CC-CCH	0	43	2	0	26	1
	CT-Kidney	0	27	1	0	31	2
AnoDet	COVID19	0	29	2	0	28	2
	INSTANCE2022	0	34	3	0	38	3
	MSD Pancreas	0	-	-	0	27	2
	ChestX-Det10	0	31	3	0	30	2
RepGene	CT-RATE	0	29	2	0	-	-
	BrainGnome	0	38	4	0	-	-
	IU_Xray	0	31	1	0	30	2
OrgSeg	BTCV	0	29	1	0	30	3
	VerSe	0	34	1	0	-	-
	L2R-OASIS	0	-	-	0	35	2

**Table 11.** Performance metrics of the Execution without collaboration (single agent) across various medical imaging tasks.

w/o Colab		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	34	4	0	-	-
	KneeMRI	0	37	5	0	-	-
	CC-CCH	0	-	-	0	-	-
	CT-Kidney	0	-	-	0	29	4
AnoDet	COVID19	0	40	5	0	-	-
	INSTANCE2022	0	36	4	0	-	-
	MSD Pancreas	0	-	-	0	33	5
	ChestX-Det10	0	41	5	0	-	-
RepGene	CT-RATE	0	-	-	0	-	-
	BrainGnome	0	32	5	0	-	-
	IU_Xray	0	-	-	0	35	4
OrgSeg	BTCV	0	-	-	0	33	5
	VerSe	0	36	5	0	-	-
	L2R-OASIS	0	32	4	0	-	-

**Table 12.** Performance metrics of the Execution without extensive debugging across various medical imaging tasks.

w/o debug		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	-	-	0	-	-
	KneeMRI	0	-	-	1	33	1
	CC-CCII	0	-	-	0	-	-
	CT-Kidney	0	-	-	0	-	-
AnoDet	COVID19	0	-	-	0	-	-
	INSTANCE2022	1	22	1	0	-	-
	MSD Pancreas	0	-	-	1	26	1
	ChestX-Det10	1	24	1	0	-	-
RepGene	CT-RATE	0	-	-	0	-	-
	BrainGnome	0	-	-	0	-	-
	IU_Xray	0	-	-	0	-	-
OrgSeg	BTCV	0	-	-	1	31	1
	VerSe	1	30	1	0	-	-
	L2R-OASIS	0	-	-	0	-	-

**Table 13.** Performance metrics of the Execution without reflection mechanism across various medical imaging tasks.

w/o reflect		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	1	37	4	0	-	-
	KneeMRI	1	35	3	1	34	5
	CC-CCII	1	31	4	1	36	4
	CT-Kidney	1	36	3	1	39	4
AnoDet	COVID19	1	25	5	1	23	5
	INSTANCE2022	1	20	3	1	23	4
	MSD Pancreas	1	27	3	1	25	4
	ChestX-Det10	0	-	-	1	31	5
RepGene	CT-RATE	1	42	6	1	46	7
	BrainGnome	1	38	5	1	39	5
	IU_Xray	0	-	-	0	-	-
OrgSeg	BTCV	1	28	2	1	27	2
	VerSe	1	31	3	0	-	-
	L2R-OASIS	0	-	-	1	27	2

**Table 14.** Performance metrics of the Execution without few-shot examples across various medical imaging tasks.

w/o fewshot		Execution 1			Execution 2		
		succ	actions	debug loops	succ	actions	debug loops
DisDiag	ADNI	0	-	-	1	37	6
	KneeMRI	1	31	4	1	29	3
	CC-CCII	1	30	3	1	32	4
	CT-Kidney	0	-	-	1	35	4
AnoDet	COVID19	1	21	3	0	-	-
	INSTANCE2022	0	-	-	0	-	-
	MSD Pancreas	1	23	3	1	26	3
	ChestX-Det10	0	-	-	1	24	4
RepGene	CT-RATE	0	-	-	0	-	-
	BrainGnome	1	33	2	1	37	4
	IU_Xray	0	-	-	1	32	3
OrgSeg	BTCV	0	-	-	1	39	5
	VerSe	1	36	4	0	-	-
	L2R-OASIS	1	37	3	0	-	-

**Table 15.** Performance metrics for different roles across tasks in OrganSeg, including Task Manager, Data Engineer, Module Architect, and Model Trainer.

OrgSeg	Task Manager				Data Engineer				Module Architect				Model Trainer			
	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn
1	1(1)	3	1	4k	1(1)	8	1	49k	1(1)	8	1	37k	1(1)	10	3	139k
2	1(1)	2	1	4k	1(1)	12	1	85k	1(1)	7	1	35k	1(1)	9	2	75k
3	1(1)	2	1	5k	0(1)	9	1	55k	1(1)	12	3	88k	0(1)	11	5	113k
4	1(1)	2	1	4k	1(1)	11	2	65k	0(1)	10	2	60k	0(1)	12	5	154k
5	1(1)	2	1	5k	1(1)	11	1	96k	1(1)	10	2	62k	1(1)	13	3	132k
6	1(1)	2	1	4k	1(1)	10	2	82k	1(1)	16	5	164k	1(1)	9	2	78k

**Table 16.** Performance metrics for different roles across tasks in AnoDet, including Task Manager, Data Engineer, Module Architect, and Model Trainer.

AnoDet	Task Manager				Data Engineer				Module Architect				Model Trainer			
	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn
1	1(1)	2	1	4k	1(1)	11	1	71k	1(1)	9	1	43k	1(1)	7	1	51k
2	1(1)	2	1	4k	1(1)	10	1	60k	1(1)	8	1	39k	1(1)	8	2	60k
3	1(1)	2	1	5k	1(1)	11	2	131k	1(1)	11	2	64k	1(1)	10	3	73k
4	1(1)	2	1	4k	1(1)	10	1	128k	0(1)	17	5	103k	1(1)	9	2	68k
5	1(1)	2	1	5k	1(1)	11	2	71k	1(1)	10	2	69k	0(1)	13	5	92k
6	1(1)	2	1	4k	1(1)	10	2	70k	1(1)	10	3	72k	1(1)	9	1	72k
7	1(1)	2	1	4k	1(1)	8	1	107k	1(1)	9	1	45k	1(1)	8	1	59k
8	1(1)	2	1	5k	1(1)	9	1	99k	1(1)	11	2	50k	1(1)	10	2	62k

**Table 17.** Performance metrics for different roles across tasks in DisDiag, including Task Manager, Data Engineer, Module Architect, and Model Trainer.

DisDiag	Task Manager				Data Engineer				Module Architect				Model Trainer			
	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn
1	1(1)	2	1	5k	1(1)	7	1	153k	1(1)	9	1	39k	1(1)	30	4	754k
2	1(1)	2	1	4k	1(1)	8	1	103k	1(1)	10	2	44k	1(1)	13	2	95k
3	1(1)	2	1	4k	1(1)	11	2	140k	1(1)	11	2	74k	1(1)	19	4	330k
4	1(1)	2	1	4k	1(1)	9	1	112k	1(1)	10	2	69k	1(1)	9	1	100k
5	1(1)	2	1	5k	1(1)	10	2	161k	0(1)	23	5	255k	1(1)	8	2	79k
6	1(1)	2	1	4k	1(1)	8	1	92k	1(1)	8	1	66k	1(1)	8	1	79k
7	1(1)	2	1	4k	1(1)	8	2	89k	1(1)	10	2	81k	1(1)	7	1	66k
8	1(1)	2	1	5k	1(1)	10	1	77k	1(1)	9	1	70k	1(1)	10	1	82k

**Table 18.** Performance metrics for different roles across tasks in RepGene, including Task Analyzer, Data Engineer, Code Writer, and Model Trainer.

RepGene	Task Analyzer				Data Engineer				Code Writer				Model Trainer			
	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn	Run	Act	Iter	Tkn
1	1(1)	2	1	4k	1(1)	4	1	14k	1(1)	5	1	18k	0(1)	14	5	103k
2	1(1)	2	1	5k	1(1)	7	1	22k	1(1)	10	2	33k	1(1)	10	2	41k
3	1(1)	2	1	4k	1(1)	11	1	94k	1(1)	16	5	171k	1(1)	9	2	80k
4	1(1)	2	1	4k	1(1)	10	1	80k	1(1)	8	1	47k	1(1)	11	1	62k
5	1(1)	2	1	4k	1(1)	11	3	110k	0(1)	17	5	219k	1(1)	13	2	79k
6	1(1)	2	1	4k	1(1)	9	1	77k	1(1)	9	1	58k	1(1)	10	1	55k