
HELIXVS: DEEP LEARNING ENHANCED STRUCTURE-BASED VIRTUAL SCREENING PLATFORM FOR HIT DISCOVERY

TECHNICAL REPORT

Shanzhuo Zhang, Xianbin Ye, Donglong He, Yueyang Huang, Xiaonan Zhang, Xiaomin Fang*
PaddleHelix Team, Baidu Inc.

August 15, 2025

ABSTRACT

Drug discovery through virtual screening (VS) has become a popular strategy for identifying hits against protein targets. VS has the potential to reduce the cost and time needed for manual selection and wet-laboratory experiments, thereby accelerating drug discovery pipelines. Improving the cost-effectiveness of virtual screening is a significant challenge, aiming to explore larger compound libraries while maintaining lower screening costs. Here, we present HelixVS, a structure-based VS platform enhanced by deep learning models. HelixVS integrates a precise deep learning-based pose-scoring model and a pose-screening module into a multi-stage VS process, enabling more effective screening of active compounds. Compared to classic molecular docking tools like Vina, HelixVS demonstrated significantly improved screening performance across nearly a hundred targets, achieving an average 2.6-fold higher enrichment factor (EF) and more than 10 times faster screening speed. We applied HelixVS in four drug development pipelines, targeting both traditional competitive drug-binding pockets and novel protein-protein interaction interfaces. Wet-lab validations across these pipelines consistently identified active compounds, with over 10% of the molecules tested in wet labs demonstrating activity at μM or even nM levels. This demonstrates the ability of HelixVS to identify high-affinity ligands for various targets and pockets. Furthermore, we provide a publicly available and free version of HelixVS with limited computing power to assist drug development scientists in accelerating their drug discovery processes. The HelixVS online service is available at: <https://paddlehelix.baidu.com/app/drug/helixvs/forecast>. We also support private deployment solutions to meet the data security requirements of pharmaceutical companies, research institutes, and other organizations. For collaboration inquiries, please contact: baidubio_cooperate@baidu.com.

Keywords Virtual screening · Deep learning · Structure-Based Drug Discovery

1 Introduction

Drug discovery is a challenging and complex process that requires extensive knowledge. High-throughput screening (HTS) is a laboratory technique used in drug discovery to identify compounds that show activity against the target of interest. Compared with HTS, a computational approach, virtual screening (VS) vastly accelerates early-stage drug discovery and development due to its cost-effectiveness and time efficiency. Virtual screening is able to explore a large chemical space containing millions to billions of compounds.

A virtual screening pipeline typically involves three main steps: (1) preparation of the target protein and compound library, (2) molecular docking and scoring to rank candidate compounds, and (3) compound filtering based on expert-defined rules. Among these, step (2) is the key step of the drug discovery pipeline. Conventional molecular docking tools, such as AutoDock, AutoDock Vina Trott and Olson [2010], Eberhardt et al. [2021], and Glide Friesner et al. [2004], Halgren et al. [2004], are used to predict ligand binding modes and affinities to target proteins. These tools are founded on the physical principles of molecular interactions, using empirical and/or physics-based scoring functions to

*Corresponding author. Email: fangxiaomin01@baidu.com

evaluate the binding affinities, and can effectively narrow down the number of potential drug candidates and prioritize those with the highest probability of success. However, they are limited by the accuracy of their scoring functions, which can result in wrong docking poses and faulty active or decoy decisions.

On the other hand, advanced deep learning-based drug-target affinity (DTA) prediction models, especially structure-based models Li et al. [2021] Jiménez et al. [2018] Jiang et al. [2021] Moon et al. [2022] Shen et al. [2022], can provide more accurate affinity predictions. They learn from large datasets of known target-ligand interactions to provide accurate predictions of binding affinities. However, the robustness of DTA model performance across different protein targets cannot be guaranteed, as the trained models may suffer from overfitting, where the model performs well on training data but fails to generalize to new data. Since both docking tools and DTA models have their strengths and weaknesses, it is attractive to integrate these two technologies into different stages of the virtual screening pipeline to take advantage of their respective benefits and improve screening accuracy.

Based on this understanding, we developed HelixVS, a cost-effective virtual screening platform to enhance screening efficiency and success rate in early-stage drug discovery. HelixVS performs multi-stage screening, comprehensively taking advantage of both molecular docking tools and deep learning-based models to screen out potential active molecules against target proteins. HelixVS demonstrated significantly improved screening performance compared to using only molecular docking tools. Compared to Vina across over a hundred targets in the DUD-E dataset Mysinger et al. [2012], HelixVS can find an average of 159% more active molecules and run up to nearly 15 times faster. These results highlight the effectiveness of HelixVS's multi-stage screening process.

Our application of the HelixVS platform in four real drug development pipelines exemplifies this approach. The targets, ranging from cyclin-dependent kinases CDK4/6 and NF- κ B inducing kinase (NIK) to immune modulators TLR4/MD-2 and cGAS, represent diverse mechanisms of disease pathogenesis and potential therapeutic intervention. The use of HelixVS enabled the screening of libraries containing millions of small molecules, identifying not only potential inhibitors with significant activity at μ M or even nM concentrations but also facilitating the exploration of novel protein-protein interaction interfaces.

The HelixVS platform offers a virtual screening service for drug development scientists, providing a comprehensive virtual screening service whose key features can be summarized as follows:

- **Superior Hit Rate:** HelixVS significantly improves virtual screening hit rates through multi-stage screening strategies and integration of deep learning models. Testing on the DUD-E dataset shows that HelixVS can achieve enrichment factors (EF) of 44.205 and 26.968 at 0.1% and 1% respectively, significantly outperforming other methods. In multiple actual drug development projects, including difficult scenarios such as dual-target and protein-protein interface (PPI) binding, HelixVS has successfully identified active molecules. These results demonstrate HelixVS platform's exceptional capability in identifying active compounds.
- **Excellent Screening Throughput and Cost-effectiveness:** HelixVS leverages Baidu Cloud's CPU computing infrastructure² and CHPC's high-performance computing resources³ to achieve screening throughput exceeding 10 million molecules per day. The platform's resource dynamic scaling capabilities and multi-stage screening workflow enable cost-effective screening across different computing resource scales, with screening costs reduced to as low as 1 RMB per thousand molecules.
- **Enhanced Usability:** HelixVS offers automated protein pre-processing, large-scale built-in compound library, custom compound library, and binding mode filtering abilities. These features are accessible through a straightforward web interface, where users can input parameters easily. This functionality is particularly beneficial for medicinal chemists who are not familiar with computational tools, allowing them to complete real virtual screening tasks quickly and accurately. The screenshots of the HelixVS web interface are shown in Appendix Section 4.

2 The HelixVS Platform

2.1 Multi-staged virtual screens with HelixVS

The HelixVS platform is a comprehensive drug discovery pipeline consisting of three stages to identify potential active molecules that can bind with a specific protein target (as shown in Figure 1).

In Stage 1, HelixVS employs classical docking tools that rely on scoring functions and sampling techniques to generate drug-target binding poses. By default, we utilize AutoDock QuickVina 2 for docking Alhossary et al. [2015], as it

²<https://cloud.baidu.com>

³<https://cloud.baidu.com/product/chpc.html>

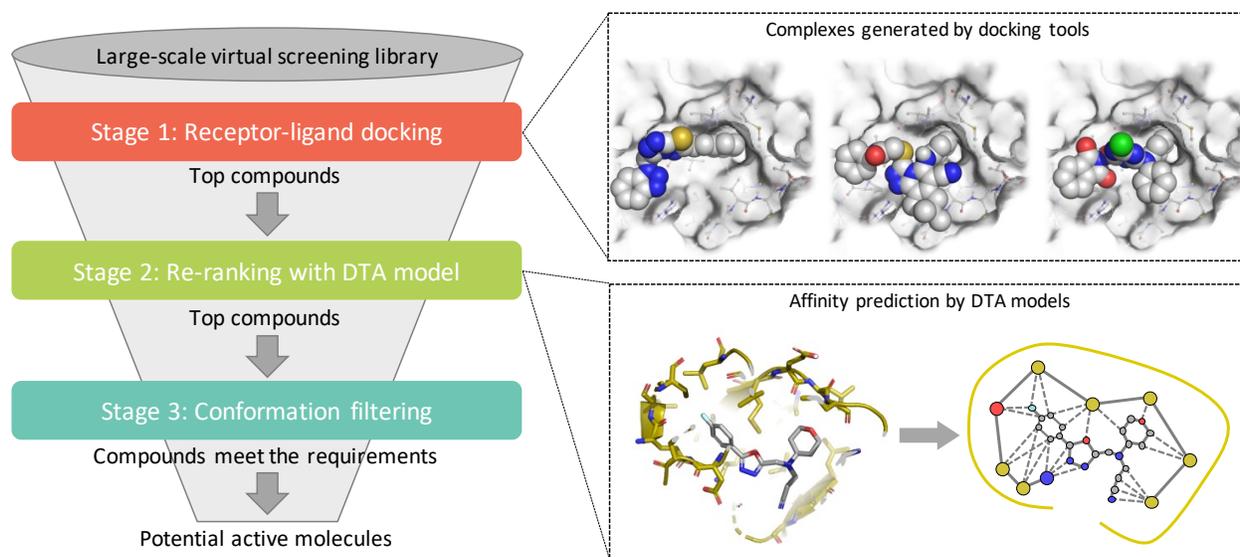


Figure 1: Overall pipeline of the HelixVS platform.

offers significantly faster speed compared to AutoDock Vina Trott and Olson [2010]. Although its scoring function is simpler and may retain non-optimal molecular conformations, we address this limitation by preserving multiple binding conformations after docking. Through our testing, we have found that retaining multiple conformations for the next stage can substantially improve enrichment factors on test datasets. Following the docking stage, molecules with higher affinity scores are selectively retained and enter the next stage.

In Stage 2, the docking poses with lower ΔG are fed into a deep learning-based affinity scoring model to obtain more accurate binding conformation scores compared to docking tools. In HelixVS, the main structure of our model is based on RTMScore, which has been proven to have high accuracy in small molecule screening scenarios Shen et al. [2022]. Building upon this foundation, we collected additional co-crystal structure data, covering almost all ligands and their spatial structures in the PDB database, to enhance RTMScore through data augmentation. At this stage, HelixVS simultaneously considers multiple isomers and uses multiple docking conformations, thereby increasing the likelihood of finding optimal conformations and affinity scores.

In Stage 3, HelixVS introduces an optional conformation filtering step to further refine molecules with higher scores based on pre-defined binding modes. This step enables targeted screening for molecules that specifically bind to certain amino acids with specific interactions. Subsequently, the remaining molecules are clustered, and representative molecules are selected to ensure diversity of the screening results.

Between each stage, HelixVS employs highly efficient distributed sorting algorithms to rank the current molecular list and gradually filter out superior molecules.

2.2 Virtual Screening Performance on DUD-E

To demonstrate the performance of HelixVS, we tested the screening power (i.e. the ability to find active molecules among a pool of active and decoy molecules) on the DUD-E Mysinger et al. [2012] dataset. DUD-E is one of the most popular virtual screening benchmarks, containing 102 proteins from 8 diverse protein families and a rich repository of 22,886 active molecules. It curates 50 topologically distinct decoys for each active molecule that possess matched physicochemical properties, sourced from the extensive ZINC database.

Here, we compare the performance of our HelixVS with other methods on DUD-E, including classic docking software Vina Trott and Olson [2010] and Glide SP Friesner et al. [2004], as well as the recently proposed deep learning-based docking model KarmaDock Zhang et al. [2023]. Implementation details are reported in Section 1.1. All methods mentioned above were evaluated on two commonly used metrics, enrichment factor Li et al. [2009] and screening speed. The overall performance is depicted in Table 1. We can see from the table that HelixVS outperforms commonly used classic docking softwares (eg. Vina and Glide SP) by a large margin on both enrichment factor and running time. Compared to Vina, HelixVS can find an average of 159% more active molecules and run up to nearly 15 times faster.

Besides, HelixVS shows an approximate 70.3% and 16.8% improvement in EF at 0.1% respectively, compared with the deep learning-based molecular docking model Karmadock and commercial software Glide SP.

Table 1: Overall performance and speed of HelixVS and other tools on DUD-E.

Method	EF _{0.1%}	EF _{1%}	Screening Speed (molecules / day)	
			per CPU core	per GPU card
Vina	17.065	10.022	~300 ^a	—
Glide SP	37.842	24.346	~2400 ^b	—
KarmaDock	25.958	15.848	—	~5 Million ^c
HelixVS	44.205	26.968	~4000	—

^a The screening speed is adapted from Gorgulla et al. [2020], with *exhaustiveness* = 32.

^b The screening speed is reported by Schrodinger [2022].

^c The screening speed is tested on NVIDIA V100 32G GPU. As reported by Zhang et al. [2023].

We also investigated the performance of all methods on each individual target families, detailed results as reported in Appendix 2.1. Surprisingly, HelixVS performs even better for certain protein categories, such as kinases, proteases, nuclear receptors, ion channels, etc. Concretely, on an important target of Parkinson’s disease, COMT (Catechol-O-methyltransferase) Ma et al. [2014], HelixVS is the only algorithm that can find active molecules, with an EF 0.1 as high as 61.51.

Recently, a new deep learning-based molecular docking model named Boltz-2 Passaro et al. [2025] has been proposed. It is a binary classification model that predicts whether a molecule can bind to a protein target based on its 3D structure. We also tested the performance of Boltz-2 on DUD-E. However, due to the significantly higher GPU resource requirements of Boltz inference compared to other methods, we selected 10 targets with smaller total ligand counts in DUD-E for inference and evaluation (specific selected targets and molecule counts are detailed in Appendix Table 2.9). The results are reported in Table 2. Boltz-2 has two different scoring modes: in binary mode, its EF 0.1 is lower than HelixVS, while EF 1 is higher than HelixVS; in affinity mode, both EF 0.1 and EF 1 are significantly lower than HelixVS. When considering computational costs, given that the cost of Nvidia GPUs with over 40GB VRAM can be 100-200 times that of a single CPU core, HelixVS’s total cost is estimated to be only 1/500 of Boltz-2’s, demonstrating HelixVS’s high cost-effectiveness in virtual screening scenarios.

Table 2: Overall performance and speed of HelixVS, Boltz-2 and other tools on 10 targets of DUD-E.

Method	EF _{0.1%}	EF _{1%}	Screening Speed (molecules / day)	
			per CPU core	per GPU card
Vina	10.948	9.308	~300	—
Glide SP	24.492	21.387	~2400	—
KarmaDock	8.924	4.104	—	~5 Million
Boltz-2 binary	33.767	27.753	—	~1440 ^a
Boltz-2 affinity	19.500	14.281	—	~1440 ^a
HelixVS	35.306	24.413	~4000	—

^a The screening speed is tested on NVIDIA A800 80G GPU. This test measures the inference speed excluding MSA retrieval, as MSA results can be reused in screening scenarios. The Boltz-2 speed is calculated based on approximately 1 minute required for inference per molecule.

To summarize, all of the above results demonstrate that HelixVS is indeed an effective virtual screening tool compared to classical docking software and deep learning-based docking model.

2.3 Compound Library Pre-processing

HelixVS supports virtual screening with both built-in compound libraries and user-uploaded custom molecular libraries, providing flexibility for diverse drug discovery needs. The built-in libraries consist of four different scales of small molecule libraries, ranging from tens of thousands to tens of millions of compounds, sourced from TopScience⁴, a leading provider of chemical compounds. This ensures that the majority of molecules in our built-in libraries are commercially available and purchasable for follow-up experimental validation. The screening time and price of built-in libraries are reported in Appendix Table 3.1. For custom molecular libraries, users can upload their own compound collections in SMILES format, enabling personalized virtual screening campaigns tailored to specific research requirements. To ensure structural rationality and computational efficiency, all compounds (both built-in and custom) undergo rigorous preprocessing steps including: (1) desalting and neutralization to standardize molecular forms, (2) stereoisomer and tautomer enumeration to account for all possible molecular conformations, (3) protonation to reflect physiological conditions, and (4) structure generation and energy filtering to eliminate energetically unfavorable conformations. This comprehensive preprocessing pipeline ensures that only high-quality, drug-like molecules with reasonable structural properties are included in the screening process, thereby improving the overall success rate of virtual screening campaigns.

3 HelixVS in Real Drug Discovering Pipelines

3.1 Screening of CDK4/6 Dual-target Inhibitors with HelixVS

Cyclin-dependent kinases 4 and 6 (CDK4/6) are key regulators of the cell cycle, controlling the progression of cells from the G1 phase to the S phase where DNA replication occurs. CDK4/6 activity is tightly regulated by the binding of cyclin D1 (CCND1), which activate the kinase domain and initiate downstream signaling pathways O'Leary et al. [2016]. Dysregulation of CDK4/6 has been implicated in the development and progression of many human cancers, making it an important target for cancer therapy. However, existing CDK4/6 inhibitors have limited specificity, leading to off-target effects and toxicity Liu et al. [2025]. In addition, cancer cells can develop resistance to these inhibitors through mutations in the ATP binding pocket or upregulation of alternative signaling pathways.

Targeting the CDK4/6-CCND1 interaction represents an alternative and more favorable strategy for inhibiting CDK4/6 activity. By preventing the formation of active CDK4/6-CCND1 complexes, these inhibitors can selectively target cancer cells with high levels of CCND1 expression, while sparing normal cells with low CCND1 expression. In addition, these inhibitors may be less susceptible to resistance mechanisms involving mutations in the ATP binding pocket Morin [2000].

Aiming at discovering a lead compound targeting CDK4/6-CCND1 interaction, HelixVS platform was used to screen a 7.8 million compound library for potential molecules that could bind to a previously unreported pocket on the CDK4/6-CCND1 interface. From the final set of 100 compounds, we purchased 40 commercially available compounds for further activity evaluation. We used an unpublished method named the bimolecular fluorescence complementation (BiFC) assay to confirm the activity of these compounds against CDK4/6-CCND1 interaction (details in the Methods section). The BiFC assay identified 6 compounds disrupting the formation of CDK4/6-CCND1 complex, as evidenced by a more than 20% decrease in fluorescence intensity. Our results demonstrate the accuracy and powerfulness of HelixVS in hit identification, underscores the platform's ability to tackle complicated cases such as the discovery of PPI inhibitors. We believe HelixVS has a great potential to make significant contribution to the development of novel therapeutics and ultimately benefit patients in need.

3.2 Screening of TLR4/MD-2 Inhibitors with HelixVS

Toll-like receptor 4 (TLR4) from the Toll-like receptor family is an important pattern recognition receptor that plays a key role in the immune system, especially in innate immune responses. It is physically associated with MD-2 and can recognize the major component of bacterial cell walls, LPS, activate downstream signaling pathways, lead to the production of inflammatory factors, and initiate inflammatory responses Shimazu et al. [1999]. Research suggests that inhibiting the activation of TLR4 may help reduce inflammatory responses in the treatment of autoimmune diseases such as rheumatoid arthritis and systemic lupus erythematosus Liu et al. [2014]. In addition, the overactivation of TLR4 also appears in some infectious diseases, allergic diseases, and cancers, making TLR4 a potential target for treatment strategies of many diseases Giat et al. [2017].

However, as of now, research on TLR4 as a drug target is mainly focused on the preclinical and clinical trial phases, and many of them are macromolecular drugs. We are committed to developing new small molecule inhibitors targeting

⁴<https://www.targetmol.cn/topscience-database>

the innovative pocket on the PPI surface between TLR4 and MD-2, which are more convenient for oral administration. Utilizing HelixVS, we meticulously screened a library of 200k molecules and selected 103 candidates for wet-lab experiments. These molecules were assessed for activity using the SEAP assay Slivka et al. [2009], which led to the identification of 6 highly potential candidates, 2 of which demonstrated nM-level activity. Follow-up modification and molecular dynamics validation of one of these molecules have been published Jiang et al. [2024].

3.3 Screening of cGAS Inhibitors with HelixVS

The cyclic GMP-AMP synthase (cGAS) acts as a cytoplasmic enzyme activated by viral infections or DNA damage, triggering downstream immune responses. Consequently, inhibiting cGAS or the STING pathway is being explored as a therapeutic strategy for autoimmune diseases like systemic lupus erythematosus and rheumatoid arthritis Hu et al. [2022]. In this project, targeting the ATP-binding pocket of cGAS, we utilized HelixVS to screen a smaller library of 30,000 compounds and identified 126 candidates. A cell-based luciferase assay Lama et al. [2019] validated 96 molecules, revealing 17 with significant activity. Remarkably, 10 of these molecules demonstrated activity below 10 μM , and one demonstrated nM-level activity.

3.4 Screening of NIK Inhibitors with HelixVS

NF- κB inducing kinase (NIK), also known as MAP3K14, is a crucial regulatory kinase in the non-canonical NF- κB signaling pathway. Clinical observations have linked the overactivation of NIK to the pathogenesis of inflammatory diseases, B-cell malignancies, and solid tumors, positioning NIK inhibition as an appealing drug discovery strategy. Such inhibitors could potentially treat conditions like cancer, inflammatory disorders, metabolic dysfunctions, and autoimmune diseases Pflug and Sitcheran [2020]. In this project, targeting the ATP-binding pocket of NIK, we utilized HelixVS to conduct a virtual screening of approximately 10 million compounds from a consolidated library to discover new active molecular scaffolds. The top-ranked molecules underwent molecular structure-based clustering, from which 7 candidates were selected for enzymatic activity assays, revealing one molecule with an IC_{50} in the μM range.

4 Conclusion

The integration of HelixVS into our drug discovery pipeline represents a significant advancement in the quest for new therapeutic agents. By utilizing this sophisticated virtual screening platform, we have successfully identified potential inhibitors that target a variety of critical proteins involved in disease mechanisms, such as CDK4/6-CCND1, TLR4/MD-2, cGAS, and NIK. These findings highlight the platform's robust capability to navigate complex molecular landscapes and pinpoint candidates with promising biological activity, some exhibiting potent inhibitory effects at μM or even nM concentrations.

Our studies have demonstrated that HelixVS not only enhances the efficiency of the screening process but also significantly improves the success rate by leveraging both classical molecular docking tools and advanced deep learning-based models. This dual approach ensures a comprehensive exploration of chemical space and a higher likelihood of discovering viable hit compounds. Furthermore, the platform's ability to handle vast libraries and deliver rapid screening results has been instrumental in accelerating the early stages of drug development.

Appendices

1 Detailed Experimental Settings

1.1 Screening Powers Evaluation

The docking protocol for Vina is as follows: the docking box is centered at the geometric center of the native ligand with a 15 Å length cubic region, followed by the standard docking pipeline, executed on CPU (Intel(R) Xeon(R) Gold 6248 CPU @ 2.50GHz) with its default scoring function and parameters except that the *exhaustiveness* is set to 32. For Glide SP, all results including docking poses and scores are freely accessible through Shen et al. [2022]. For KarmaDock, we load the trained model parameters provided by KarmaDock Zhang et al. [2023] and follow their settings to report the results, using NVIDIA V100 32G GPU.

1.2 The Bimolecular Fluorescence Complementation Assay

The BiFc assay was used to detect protein-protein interactions between CDK4/6 and CCND1. Briefly, the sequences of VN155-CDK4/6 and VC155-CCND1 were cloned into the vectors carrying the N- and C-terminal halves of the green fluorescent protein (GFP), respectively. The constructed vectors were confirmed by sequencing. HEK293T cells were grown to 50% confluency and co-transfected with VN155-CDK4/VC155-CCND1 or VN155-CDK6/VC155-CCND1 vectors. The PPI fluorescence signal was detected by fluorescence microscopy after 24 hours (Figure). To evaluate the inhibitory activities of the tested compounds, transfected cells were treated with compounds at a final concentration of 10 μM for 18 hours. The fluorescence intensity was quantified using a high-content imaging system. Cells treated with DMSO was used as a blank. The mCherry signal intensity was used as a reference signal for normalizing the PPI signal.

2 Detailed Results on DUD-E

2.1 Virtual Screening Performance about each individual target family on DUD-E

We report the performance (EF 0.1% and EF 1%) of each method on an individual target family, as shown from table 2.1 to 2.8. In short, HelixVS achieved state-of-the-art results on EF 0.1% on six of the eight target families, compared to the other methods. Especially in the Ion Channel target family, HelixVS surpasses all other baselines by a large margin, achieving a 5.847-fold and 2.560-fold improvement compared with Vina and Glide SP, respectively. Unfortunately, due to the complexity of the binding mechanism of the GPCR and Cytochrome P450 targets, HelixVS performed slightly less effectively than the best methods, finishing as the runner-up. Considering its robust performance, we can expect the promising prospect of HelixVS in real-world virtual screening tasks.

Table 2.1: Cytochrome P450.

Method	EF _{0.1%}	EF _{1%}
Vina	3.256	3.607
Glide SP	6.399	3.923
KarmaDock	12.732	4.347
HelixVS	6.512	4.811

Table 2.2: GPCR.

Method	EF _{0.1%}	EF _{1%}
Vina	3.966	2.186
Glide SP	31.555	19.129
KarmaDock	6.298	1.736
HelixVS	23.190	13.616

Table 2.3: Ion Channel.

Method	EF _{0.1%}	EF _{1%}
Vina	10.348	7.734
Glide SP	19.905	14.200
KarmaDock	11.007	1.122
HelixVS	70.854	42.150

Table 2.4: Kinase.

Method	EF _{0.1%}	EF _{1%}
Vina	22.453	13.118
Glide SP	48.765	28.640
KarmaDock	45.624	30.565
HelixVS	51.890	30.183

Table 2.5: Miscellaneous.

Method	EF _{0.1%}	EF _{1%}
Vina	4.002	10.887
Glide SP	33.325	26.342
KarmaDock	8.971	3.735
HelixVS	41.325	23.661

Table 2.6: Nuclear Receptor.

Method	EF _{0.1%}	EF _{1%}
Vina	27.834	13.672
Glide SP	48.071	29.355
KarmaDock	16.227	8.752
HelixVS	48.518	25.617

Table 2.7: Other Enzymes.

Method	EF _{0.1%}	EF _{1%}
Vina	16.374	8.933
Glide SP	27.571	20.084
KarmaDock	13.459	8.056
HelixVS	36.031	25.300

Table 2.8: Protease.

Method	EF _{0.1%}	EF _{1%}
Vina	13.307	8.081
Glide SP	46.242	28.609
KarmaDock	44.980	26.479
HelixVS	57.286	32.871

Table 2.9: 10 targets selected from DUD-E for Boltz-2 evaluation

Target Name	Class	#Actives	#Decoys
AMPC	Other Enzymes	48	2850
COMT	Other Enzymes	41	3850
GLCM	Other Enzymes	54	3800
INHA	Other Enzymes	43	2300
KITH	Other Enzymes	57	2850
PUR2	Other Enzymes	50	2700
PYGM	Other Enzymes	77	3950
SAHH	Other Enzymes	63	3450
FABP4	Miscellaneous	47	2750
CXCR4	GPCR	40	3406

3 Built-in compound libraries and screening pricing

4 Screenshots of HelixVS Server

Table 3.1: Built-in compound libraries and screening pricing

Library Name	#Molecules	Screening Time (h) ^a	Screening Price (RMB)
Targetmol_CherryPick	~34,000	<1	70
Lifechemicals	~820,000	<6	1,200
ChemDiv	~2,830,000	<11	3,700
Topscience database	~16,870,000	<26	19,000

^a Screening time when resources are uncontested, may occur when multiple "Topscience database" tasks are running.

HelixVS

- High cost-effectiveness virtual screening process
- Supports billion-level virtual screening (Form filling required)
- Hit rate increased by up to 159%
- Up to 28 times acceleration

[Tutorial](#)

1 Protein preparation 2 Define binding pocket 3 Set screening library and parameters


Data input Search in database Upload a file Use an example

Fill with an example 2FSZ

```

Sequence of 2FSZ.pdb Chain 1: Polymer 1 (... A)
264 274 284 294 304 314 324
SPFQVLTLLLEAPPVHLISRPGAPFTEASMMMSLTKLADRELVHMISWAKIIPGFVLSLFDQV
324 344 354 364 374 384 39
RLESSMVEVLMGMLMWSIDHPGKLI FAFDVLDRDEGKSVGILEIFDMLLATTSRFRELKIQ
404 426 436 446 456 466
HKEYLCVKAMILLNSSMSRKLHLNVAVDALVWVIKSGISSQQSMRLAMLLMLLSHVRHASN

```

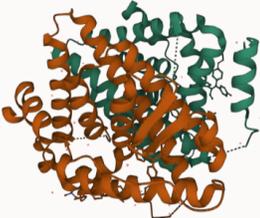

Next

Figure 4.1: The protein preparation page. The protein structure can be uploaded in PDB format or downloaded from PDB database.

HelixVS

- High cost-effectiveness virtual screening process
- Supports billion-level virtual screening (Form filling required)
- Hit rate increased by up to 159%
- Up to 28 times acceleration

[Tutorial](#)

1 **Protein preparation** 2 **Define binding pocket** 3 **Set screening library and parameters**

Target protein: 2FSZ

Data input:

```
Sequence of 2FSZ.pdb Chain 1: Polymer 1 (... A
264 274 284 294 304 314 324
SPEQLVLTLEAEPFHVLI SRPSAPFTEASMMMSLTKLADKELVHMISWAKI PGFVELSLFDQV
324 334 344 354 364 374 384 39
RLESSEVLMGMLMRSIDHPGKLIFADPDLVDRDEGKSVGEILEIFDMLLATTSRFRELKIQ
404 414 424 434 444 454 464
HKEYLCVKAMILLNSSMRKLAHLNNAVTDALVWVIKSGISQQSMRLANLMLLSHVREASN
```

Center coordinates
X: Y: Z:

Pocket size
W: H: L:

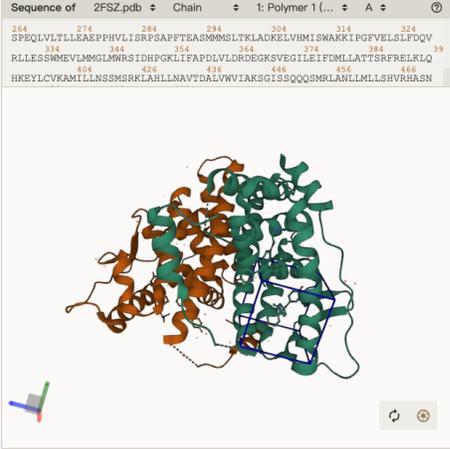


Figure 4.2: The pocket specification page. The pocket can be specified by setting a box with center coordinates and size on three dimensions in Angstrom. The pocket will be displayed on the protein structure by clicking the "Show Pocket" button.

HelixVS

- High cost-effectiveness virtual screening process
- Supports billion-level virtual screening (Form filling required)
- Hit rate increased by up to 159%
- Up to 28 times acceleration

[Tutorial](#)

1 Protein preparation 2 Define binding pocket **3 Set screening library and parameters**

Data input:

Select screening library: Targetmol_CherryPick (34,022)

Lifechemicals (820,215)

ChemDiv (2,830,000)

topscience database (16,870,000)

Maximum output molecule count:

Interaction constraints:

Task name:

Figure 4.3: The parameter setting page. The user can set the screening library, and the interaction constraints.

Results for HelixVS Task name: demo Task ID:103198

Get free credits by evaluation ★★★★★

Screening library: Targetmol_CherryPick (34,022)

Maximum output molecule count: 25

Interaction constraints: -

Molecule	DTI ...	Sequence of 2X05	Chain
	53.815	AGNATP1FGVVIQSSAQVSDDAVSKPGFPTSGWYFVSSRSTVYAGLLQNGKYADPFYSTNMQNVFAAQFVPMWYRDLAVDQTSRRTYLDVDFSGVLSKADVVNGTKVATKQVNGAYTBHDLDTAQVHTGVNSVAFKYPNDPNRDLMSGWIW	1: Polymer 1
	50.451	AQTFDQNGMIVEDLVRRSGVALRSRHHVIQKLNLSALDHADLVKADVANDSANAQVTVAGTVAGKPISTQVSLAAKERKTVTFPLVGLDRFNWPPAGMGQQRHYLDLDTASVGGTSDAAKSKFGVRDVKATLMSGGQRYVNGKPLLIQGG	1: Polymer 1
	50.069	GYTFDLFRWNETAAAKLKYVNLGLNIVRLEGHIPEDEFFIADDGLVLTMPGWGCCDRWEGQVNGEKGEPWVESDYPIAKASMFSEARLRDHPVSI SFHIGSDFAPDRRIEQGYLDAMKAADFLLPVIPAASARFSPITGASGMQNNNGPYDY	1: Polymer 1
	48.926		
	48.778		

Show pocket

Figure 4.4: The result page showing the screening task results. The left panel displays the filtered molecular list and 2D topological molecular structures, while the right panel shows the binding structure of the currently selected molecule with the target protein. Structure visualization is implemented using Mol* Viewer Sehnal et al. [2021].

References

- Oleg Trott and Arthur J Olson. Autodock vina: improving the speed and accuracy of docking with a new scoring function, efficient optimization, and multithreading. *Journal of computational chemistry*, 31(2):455–461, 2010.
- Jerome Eberhardt, Diogo Santos-Martins, Andreas F. Tillack, and Stefano Forli. Autodock vina 1.2.0: New docking methods, expanded force field, and python bindings. *Journal of Chemical Information and Modeling*, 61(8):3891–3898, 2021. doi:10.1021/acs.jcim.1c00203. URL <https://doi.org/10.1021/acs.jcim.1c00203>. PMID: 34278794.
- Richard A Friesner, Jay L Banks, Robert B Murphy, Thomas A Halgren, Jasna J Klicic, Daniel T Mainz, Matthew P Repasky, Eric H Knoll, Mee Shelley, Jason K Perry, et al. Glide: a new approach for rapid, accurate docking and scoring. 1. method and assessment of docking accuracy. *Journal of medicinal chemistry*, 47(7):1739–1749, 2004.
- Thomas A Halgren, Robert B Murphy, Richard A Friesner, Hege S Beard, Leah L Frye, W Thomas Pollard, and Jay L Banks. Glide: a new approach for rapid, accurate docking and scoring. 2. enrichment factors in database screening. *Journal of medicinal chemistry*, 47(7):1750–1759, 2004.
- Shuangli Li, Jingbo Zhou, Tong Xu, Liang Huang, Fan Wang, Haoyi Xiong, Weili Huang, Dejing Dou, and Hui Xiong. Structure-aware interactive graph neural networks for the prediction of protein-ligand binding affinity. In *Proceedings of the 27th ACM SIGKDD Conference on Knowledge Discovery & Data Mining*, pages 975–985, 2021.
- José Jiménez, Miha Skalic, Gerard Martinez-Rosell, and Gianni De Fabritiis. K deep: protein–ligand absolute binding affinity prediction via 3d-convolutional neural networks. *Journal of chemical information and modeling*, 58(2): 287–296, 2018.
- Dejun Jiang, Chang-Yu Hsieh, Zhenxing Wu, Yu Kang, Jake Wang, Ercheng Wang, Ben Liao, Chao Shen, Lei Xu, Jian Wu, et al. Interactiongraphnet: A novel and efficient deep graph representation learning framework for accurate protein–ligand interaction predictions. *Journal of medicinal chemistry*, 64(24):18209–18232, 2021.
- Seokhyun Moon, Wonho Zhung, Soojung Yang, Jaechang Lim, and Woo Youn Kim. Pignet: a physics-informed deep learning model toward generalized drug–target interaction predictions. *Chemical Science*, 13(13):3661–3673, 2022.
- Chao Shen, Xujun Zhang, Yafeng Deng, Junbo Gao, Dong Wang, Lei Xu, Peichen Pan, Tingjun Hou, and Yu Kang. Boosting protein–ligand binding pose prediction and virtual screening based on residue–atom distance likelihood potential and graph transformer. *Journal of Medicinal Chemistry*, 65(15):10691–10706, 2022. doi:10.1021/acs.jmedchem.2c00991. URL <https://doi.org/10.1021/acs.jmedchem.2c00991>. PMID: 35917397.
- Michael M. Mysinger, Michael Carchia, John. J. Irwin, and Brian K. Shoichet. Directory of useful decoys, enhanced (dud-e): Better ligands and decoys for better benchmarking. *Journal of Medicinal Chemistry*, 55(14):6582–6594, 2012. doi:10.1021/jm300687e. URL <https://doi.org/10.1021/jm300687e>. PMID: 22716043.
- Amr Alhossary, Stephanus Daniel Handoko, Yuguang Mu, and Chee-Keong Kwoh. Fast, accurate, and reliable molecular docking with quickvina 2. *Bioinformatics*, 31(13):2214–2216, 2015.
- Xujun Zhang, Odin Zhang, Chao Shen, Wanglin Qu, Shicheng Chen, Hanqun Cao, Yu Kang, Zhe Wang, Ercheng Wang, Jintu Zhang, Yafeng Deng, Furui Liu, Tianyue Wang, Hongyan Du, Langcheng Wang, Peichen Pan, Guangyong Chen, Chang-Yu Hsieh, and Tingjun Hou. Efficient and accurate large library ligand docking with karmadock. *Nature Computational Science*, 3(9):789–804, Sep 2023. ISSN 2662-8457. doi:10.1038/s43588-023-00511-5. URL <https://doi.org/10.1038/s43588-023-00511-5>.
- Honglin Li, Hailei Zhang, Mingyue Zheng, Jie Luo, Ling Kang, Xiaofeng Liu, Xicheng Wang, and Hualiang Jiang. An effective docking strategy for virtual screening based on multi-objective optimization algorithm. *BMC Bioinformatics*, 10(1):58, Feb 2009. ISSN 1471-2105. doi:10.1186/1471-2105-10-58. URL <https://doi.org/10.1186/1471-2105-10-58>.
- Christoph Gorgulla, Andras Boeszoermyenyi, Zi-Fu Wang, Patrick D. Fischer, Paul W. Coote, Krishna M. Padmanabha Das, Yehor S. Malets, Dmytro S. Radchenko, Yurii S. Moroz, David A. Scott, Konstantin Fackeldey, Moritz Hoffmann, Iryna Iavniuk, Gerhard Wagner, and Haribabu Arthanari. An open-source drug discovery platform enables ultra-large virtual screens. *Nature*, 580(7805):663–668, Apr 2020. ISSN 1476-4687. doi:10.1038/s41586-020-2117-z. URL <https://doi.org/10.1038/s41586-020-2117-z>.
- Schrodinger. How long does it take to screen 10,000 compounds with glide ?, 2022. https://support.schrodinger.com/s/article/1012?original_search=glide%20speed, Last Published on 7/29/2022 11:47 PM.
- Zhiguo Ma, Hongming Liu, and Baojian Wu. Structure-based drug design of catechol-o-methyltransferase inhibitors for cns disorders. *British Journal of Clinical Pharmacology*, 77(3):410–420, 2014. doi:<https://doi.org/10.1111/bcp.12169>. URL <https://bpspubs.onlinelibrary.wiley.com/doi/abs/10.1111/bcp.12169>.

- Saro Passaro, Gabriele Corso, Jeremy Wohlwend, Mateo Reveiz, Stephan Thaler, Vignesh Ram Somnath, Noah Getz, Tally Portnoi, Julien Roy, Hannes Stark, David Kwabi-Addo, Dominique Beaini, Tommi Jaakkola, and Regina Barzilay. Boltz-2: Towards accurate and efficient binding affinity prediction, 2025. URL <https://www.biorxiv.org/content/10.1101/2025.06.14.659707v1>. Pages: 2025.06.14.659707 Section: New Results.
- Ben O’Leary, Richard S. Finn, and Nicholas C. Turner. Treating cancer with selective CDK4/6 inhibitors. 13(7):417–430, 2016. ISSN 1759-4782. doi:10.1038/nrclinonc.2016.26. URL <https://doi.org/10.1038/nrclinonc.2016.26>.
- Yanbiao Liu, Seohyun Park, and Yan Li. Breaking cancer’s momentum: CDK4/6 inhibitors and the promise of combination therapy. 17(12):1941, 2025. ISSN 2072-6694. doi:10.3390/cancers17121941. URL <https://www.mdpi.com/2072-6694/17/12/1941>. Number: 12 Publisher: Multidisciplinary Digital Publishing Institute.
- Michael J Morin. From oncogene to drug: development of small molecule tyrosine kinase inhibitors as anti-tumor and anti-angiogenic agents. 19(56):6574–6583, 2000. ISSN 1476-5594. doi:10.1038/sj.onc.1204102. URL <https://doi.org/10.1038/sj.onc.1204102>.
- Rintaro Shimazu, Sachiko Akashi, Hirotaka Ogata, Yoshinori Nagai, Kenji Fukudome, Kensuke Miyake, and Masao Kimoto. MD-2, a Molecule that Confers Lipopolysaccharide Responsiveness on Toll-like Receptor 4. 189(11):1777, 1999. doi:10.1084/jem.189.11.1777. URL <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC2193086/>.
- Yu Liu, Heng Yin, Ming Zhao, and Qianjin Lu. TLR2 and TLR4 in autoimmune diseases: a comprehensive review. 47(2):136–147, 2014. ISSN 1559-0267. doi:10.1007/s12016-013-8402-y. URL <https://doi.org/10.1007/s12016-013-8402-y>.
- Eitan Giat, Michael Ehrenfeld, and Yehuda Shoenfeld. Cancer and autoimmune diseases. 16(10):1049–1057, 2017. ISSN 1568-9972. doi:10.1016/j.autrev.2017.07.022. URL <https://www.sciencedirect.com/science/article/pii/S1568997217302057>.
- Peter F. Slivka, Mitesh Shridhar, Gui-in Lee, Deanne W. Sammond, Mark R. Hutchinson, Alexander J. Martinko, Madison M. Buchanan, Page W. Sholar, Jeffrey J. Kearney, Jacqueline A. Harrison, Linda R. Watkins, and Hang Yin. A peptide antagonist of the tr4–md2 interaction. 10(4):645–649, 2009. ISSN 1439-4227, 1439-7633. doi:10.1002/cbic.200800769. URL <https://chemistry-europe.onlinelibrary.wiley.com/doi/10.1002/cbic.200800769>.
- Yao-Yao Jiang, Shuai-Ting Yan, Shan-Zhuo Zhang, Meng Wang, Wei-Ming Diao, Jun Li, Xiao-Min Fang, and Hang Yin. Discovery of pyrazolo[1,5-*a*]pyrimidine derivatives targeting tr4 homodimerization via AI-powered next-generation screening. 280:116945, 2024. ISSN 0223-5234. doi:10.1016/j.ejmech.2024.116945. URL <https://www.sciencedirect.com/science/article/pii/S0223523424008262>.
- Ying Hu, Bangjie Chen, Fan Yang, Yue Su, Dashuai Yang, Yan Yao, Shuxian Wang, Yincui Wu, Liangsong Tao, and Tao Xu. Emerging role of the cGAS-STING signaling pathway in autoimmune diseases: Biologic function, mechanisms and clinical prospection. 21(9):103155, 2022. ISSN 1568-9972. doi:10.1016/j.autrev.2022.103155. URL <https://www.sciencedirect.com/science/article/pii/S1568997222001252>.
- Lodoe Lama, Carolina Adura, Wei Xie, Daisuke Tomita, Taku Kamei, Vitaly Kuryavyi, Tasos Gogakos, Joshua I. Steinberg, Michael Miller, Lavoisier Ramos-Espiritu, Yasutomi Asano, Shogo Hashizume, Jumpei Aida, Toshihiro Imaeda, Rei Okamoto, Andy J. Jennings, Mayako Michino, Takanobu Kuroita, Andrew Stamford, Pu Gao, Peter Meinke, J. Fraser Glickman, Dinshaw J. Patel, and Thomas Tuschl. Development of human cGAS-specific small-molecule inhibitors for repression of dsDNA-triggered interferon expression. 10(1):2261, 2019. ISSN 2041-1723. doi:10.1038/s41467-019-08620-4. URL <https://www.nature.com/articles/s41467-019-08620-4>.
- Kathryn M. Pflug and Raquel Sitcheran. Targeting NF- κ B-Inducing Kinase (NIK) in Immunity, Inflammation, and Cancer. 21(22):8470, 2020. ISSN 1422-0067. doi:10.3390/ijms21228470.
- David Sehnal, Sebastian Bittrich, Mandar Deshpande, Radka Svobodová, Karel Berka, Václav Bazgier, Sameer Velankar, Stephen K Burley, Jaroslav Koča, and Alexander S Rose. Mol* viewer: modern web app for 3d visualization and analysis of large biomolecular structures. 49:W431–W437, 2021. ISSN 0305-1048. doi:10.1093/nar/gkab314. URL <https://doi.org/10.1093/nar/gkab314>. _eprint: <https://academic.oup.com/nar/article-pdf/49/W1/W431/38842088/gkab314.pdf>.