

# Adaptive and robust smartphone-based step detection in multiple sclerosis

Lorenza Angelini,<sup>1</sup> Dimitar Stanev,<sup>1</sup> Marta Płonka,<sup>2</sup> Rafał Klimas,<sup>1</sup> Natan Napiórkowski,<sup>2</sup> Gabriela González Chan,<sup>3</sup> Lisa Bunn,<sup>3</sup> Paul S Glazier,<sup>3</sup> Richard Hosking,<sup>4</sup> Jenny Freeman,<sup>3</sup> Jeremy Hobart,<sup>5</sup> Jonathan Marsden,<sup>3</sup> Licinio Craveiro,<sup>1</sup> Mike D Rinderknecht,<sup>1\*</sup> and Mattia Zanon<sup>1\*</sup>

\* Shared last author

<sup>1</sup>F. Hoffmann-La Roche Ltd, Basel, Switzerland

<sup>2</sup>Roche Polska Sp. z o.o., Warsaw, Poland

<sup>3</sup>School of Health Professions, University of Plymouth, Plymouth, UK

<sup>4</sup>University of Plymouth Enterprise Ltd, Plymouth, UK

<sup>5</sup>Peninsula Medical School, Faculty of Health, University of Plymouth, Plymouth, UK

## Corresponding author

Mike D Rinderknecht

[mike.rinderknecht@roche.com](mailto:mike.rinderknecht@roche.com)

# Abstract

**Background:** Many attempts to validate gait pipelines, a framework that processes sensor data to detect gait events for subsequent computation of gait-related measures, have focused on validating the detection of initial contacts over final contacts in supervised environments using a single sensor.

**Objective:** To evaluate the performance of a gait pipeline in detecting initial and final contacts, (thus enabling computation of a larger set of gait measures) using a step detection algorithm adaptive to different testing environments, smartphone wear locations, and gait impairment levels.

**Methods:** In GaitLab, healthy controls (HC) and people with multiple sclerosis (PwMS; Expanded Disability Status Scale 0.0–6.5) performed supervised laboratory-based walking activities (structured in-lab overground and treadmill Two-Minute Walk Test [2MWT]) during two on-site visits carrying six smartphones (left/right front pocket, front/back belt and left/right back pocket) and unsupervised walking activities (structured and unstructured real-world walking) once daily for 10–14 days using a single smartphone. Reference gait data were simultaneously collected with either a motion capture system (structured in-lab treadmill 2MWT) or Gait Up inertial measurement unit sensors (all other assessments). The gait pipeline's performance in detecting gait events (initial/ final contacts) was evaluated through F1 scores (harmonic mean of precision and recall; tolerance window: 0.5 s) and absolute temporal error with respect to reference measurement systems.

**Results:** We studied 35 HC and 93 PwMS. Initial and final contacts were accurately detected across all smartphone wear locations. On in-lab walking activities, median F1 scores for initial/final contacts were  $\geq 98.2\%/96.5\%$  in HC and  $\geq 98.5\%/97.7\%$  in PwMS. Median (IQR) F1 scores remained high on structured real-world walking (HC: 100% [0.3%]/100% [0.2%]; PwMS: 99.5% [1.9%]/99.4% [2.5%]) and unstructured real-world walking (HC: 97.8% [2.6%]/97.8%

[2.8%]; PwMS: 94.4% [6.2%]/94.0% [6.5%]). Median temporal errors for both initial and final contacts were  $\leq 0.08$  s. Neither age, sex, disease severity, walking aid use, nor setting (outdoor/indoor) impacted significantly on pipeline performance (all  $p > 0.05$ ).

**Conclusion:** This gait pipeline accurately and consistently detected initial and final contacts in PwMS across different smartphone locations and environments, expanding previous research and highlighting its potential for real-world gait assessment.

**Trial registration:** ISRCTN15993728

**Keywords:** Real-world data, Gait, Validation, Multiple sclerosis, Smartphone, Two-Minute Walk Test, Gait Impairment, Sensors

## Introduction

Gait impairment is commonly observed among people with multiple sclerosis (PwMS) [1, 2]. Slowing or stopping worsening of walking disability can have a profound impact on affected patients and the burden the disease has on society [3-6]. Different scales and performance measures are available to assess and monitor gait impairment including the Two-Minute Walk Test (2MWT) [7], Timed 25-Foot Walk (T25FW) [8], Timed Up and Go (TUG) [9], and Expanded Disability Status Scale (EDSS) [10]. Despite their use, they have significant limitations in capturing an individual's functional abilities in daily life and identifying subtle changes in gait, particularly in the early stages of disease. These limitations arise from the requirement for a controlled environment and the presence of a trained examiner, which differs from the challenges that PwMS face daily [11]. Furthermore, these scales and performance measures only provide a single value, whether it is the time needed to complete the task or a score, and thus only partially reflect the complex and dynamic nature of mobility, such as the capacity to

navigate uneven surfaces or respond to unexpected obstacles [12]. Alternative assessment methods, such as wearable technologies, may offer a comprehensive and quantitative approach and have the potential to provide more detailed information on several gait-related parameters in daily life [13-19]. Smartphones, in particular, have emerged as a promising tool for frequent, cost-efficient, convenient, and minimally intrusive assessments in a home environment and provide a much larger set of measures to probe different aspects of functional ability [20-22]. Incorporating these alternative assessment methods ensures a more accurate assessment of functional impairment, enabling the development of more effective therapeutic options for PwMS.

Multiple pipelines that extract digital outcomes from gait data recorded through inertial measurement units (IMUs) embedded within smartphones or wearable devices [13, 23-26] have been developed in recent years. However, variations in sensor placement, disease characteristics, contextual, and environmental factors can significantly impact the technical validity of the step detection algorithms (i.e., the extent to which such algorithms can accurately and reliably detect steps when compared against reference measurement systems) integrated in such gait pipelines, especially in populations with pathological gait, such as PwMS. Key findings from the existing literature highlight the following issues:

- **Sensor placement:** Many validation studies have predominantly used sensors placed on a fixed location, such as the lower back or the ankles [14, 16, 27].
- **Gait events:** While some gait pipelines are available that are robust against different sensor placements, they are typically designed to count steps and are, therefore, restricted to detecting initial contacts only, thus limiting the range of digital gait measures they can compute [28, 29].
- ✶ **Population type:** Technical validation of step detection algorithms integrated in gait pipelines is mainly conducted in healthy controls (HC) [28, 30, 31], with only limited

studies addressing pathological populations, such as those with Parkinson's disease [32, 33].

- **Context and environment:** Step detection algorithms are often validated in controlled environments, such as laboratories or clinics, which do not reflect performance in everyday life [34-36].

These limitations highlight the importance of conducting validation studies that examine the potential impact of different smartphone wear locations and other contextual and environmental factors on the detection of initial and final contacts (i.e., heel strike and toe-off events, respectively) in a representative sample of individuals with relevant pathologies. Addressing these issues is critical for developing a robust gait pipeline that can reliably characterize gait abnormalities in daily life contexts across diverse patient populations, such as people with multiple sclerosis.

The GaitLab study (ISRCTN15993728) was designed to collect data for the development and validation of a gait pipeline that is robust with regards to different intrinsic gait factors (e.g., population type, disease stage, walking speeds), contextual and environment factors (e.g., indoor vs outdoor setting, in-lab vs real-world environment, structured vs unstructured walking activities, use of walking aids and orthotics), and smartphone wear locations [37]. Using data from this study collected from PwMS and HC, we establish here the technical validity of a gait pipeline and, in particular, of a novel step detection adaptive to different testing environments, smartphone wear locations, and gait impairment levels. We also describe the step detection algorithm's performance in accurately identifying the gait events within a gait cycle such as initial and final contacts. Furthermore, we evaluate the impact of intrinsic gait factors as well as contextual and environmental factors on the step detection algorithm's performance.

# Methods

## Study design and participants

The GaitLab study employed a well-defined experimental protocol and state-of-the-art equipment to develop and validate the gait pipeline. The study design has been previously described in detail [37]. In short, the study recruited PwMS and age- and sex-matched HC from a single site in Plymouth, United Kingdom, from January 2022 to October 2023. All participants were required to be at least 18 years old. PwMS were additionally required to have a confirmed diagnosis of either relapsing-remitting, primary progressive, or secondary progressive MS, per the 2010 or 2017 McDonald diagnostic criteria [38, 39]; having the ability to walk for a period of at least 6 minutes, and a baseline EDSS score between 0.0 and 6.5, inclusive. The protocol consisted of two visits to a gait laboratory and a two-week period of unsupervised remote testing in the home environment in between (Figure 1). The technical validity of the gait pipeline and in particular of the novel and adaptive step detection integrated in the pipeline was evaluated using gait data obtained during the second gait laboratory visit as well as during the remote testing period. Study participants were permitted to wear standard footwear and were allowed to use any necessary walking aids or orthotics, or rest as required (participants were asked to use the same walking aids or orthotics during all sessions if they used any).

## In-lab walking activities

During the second on-site visit, each participant was requested to perform a structured in-lab overground walking task and a structured in-lab treadmill walking task (Figure 1). For the first task, participants were instructed to complete a 2MWT, which involved walking back and forth along a marked 10-meter corridor, while walking safely at their fastest possible pace. For the

second task, participants completed the 2MWT on a dual-belt treadmill (Motek, Netherlands), which featured a self-paced algorithm to enable participants to walk at a self-selected fast pace, i.e., a pace that resembles overground walking most closely [40]. During the treadmill session, participants could hold the handrails in the event of losing balance and were outfitted with a safety harness to prevent falls.


Participants completed both activities while carrying six smartphones with a built-in IMU (Samsung Galaxy A40 running on Android version 9; sampling frequency: 50 Hz) in customized shorts and adjustable waist belts to capture gait data. Each smartphone was placed in different pockets and belts. These were the left/right front pocket, belt back, belt front, and left/right back pocket (Figure 1).

## Real-world walking activities

During the unsupervised remote testing in the home environment, participants performed unsupervised and self-administered structured and unstructured real-world walking activities in real-life conditions for 10-14 days. The structured real-world walking activity entailed a daily 2MWT, during which participants were directed to walk as fast and far as possible but safely for 2 minutes while minimizing sharp turns as much as possible, mimicking the clinician-administered 2MWT. This constraint was introduced so that this activity resembles the clinician-administered 2MWT as closely as possible. Participants were required to carry a single smartphone, which was preferably placed in a waist-worn front belt bag (Figure 1). The structured real-world walking activity could be conducted indoors or outdoors, with a preference for alternating between the two on consecutive days if possible. By comparison, the unstructured real-world walking activity involved the study participants carrying a single smartphone for up to 4 hours a day with at least 15 minutes of walking while they went about their typical daily routine.

Study participants were allowed to position the smartphone in their preferred location, ensuring it remained as close to their body as possible (Figure 1).

Following each structured and unstructured real-world walking activity, participants were required to update their daily gait diary, in which they captured contextual information about how they performed the gait tests remotely, such as the smartphone's wear location, whether walking was performed indoors or outdoors, information on the terrain (e.g., flat versus hilly) if walking outdoors, or details on walking aids if used. The gait diary was designed to explore contextual determinants of gait quality, providing a more comprehensive analysis of the collected data.

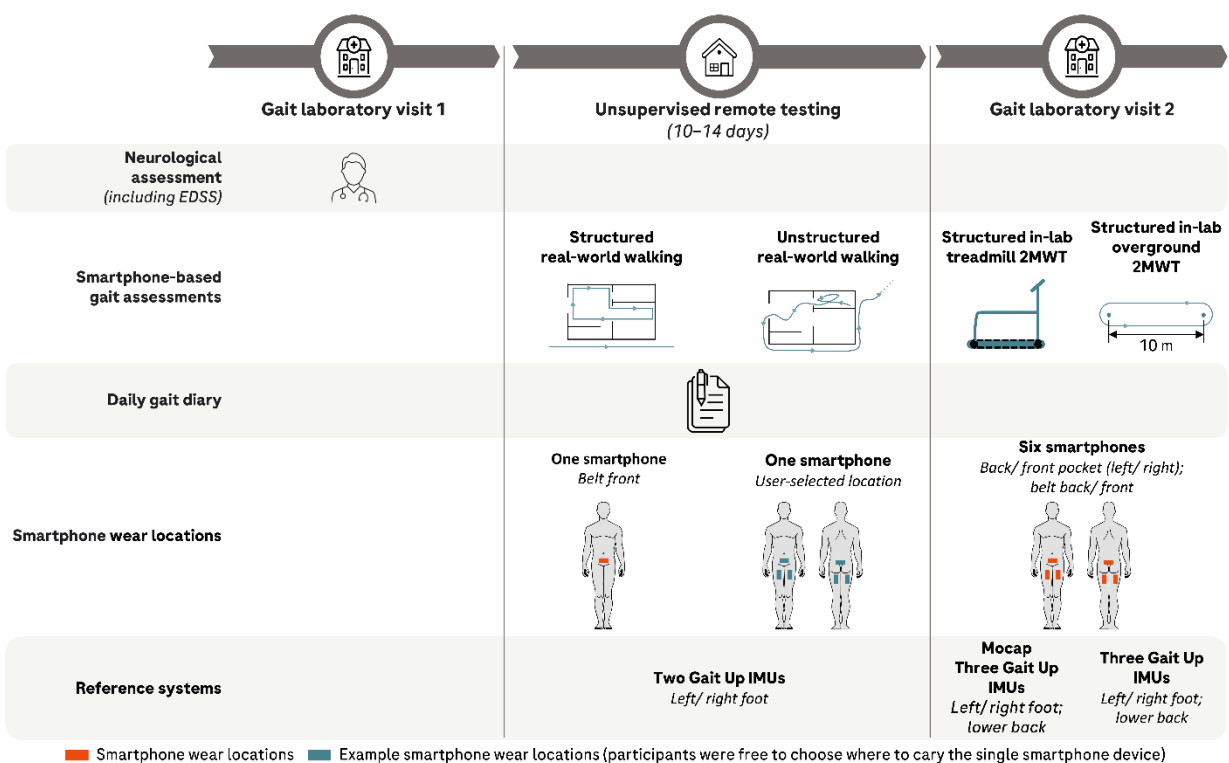

**Figure 1. Outline of the clinical visits and relevant data used in the analyses.** The gait pipeline was validated using gait data collected from the structured real-world walking, unstructured real-world walking, structured in-lab treadmill 2MWT, and structured in-lab overground 2MWT. Data collected during the neurological assessment and with the daily gait diary were used to provide additional context for the interpretation of the gait pipeline's performance metrics. A full overview of all assessments and data collected in the study is reported in [37]. 2MWT, Two-Minute Walk Test; IMUs, inertial measurement units.

## Reference measurement systems

To evaluate the step detection algorithm's technical validity, ground-truth data from two reference measurement systems were collected simultaneously with smartphone data. The reference measurement systems were synchronized to the smartphones as described in section “Supplementary appendix” and in [41]. The “gold standard” reference system consisted of a marker-based motion capture (mocap) system and two multiaxial force plates. The mocap system included 12 infrared cameras (Vicon Vero™ v1.3 cameras) with a resolution of 1.3 MP at 100 frames per second and 26 reflective markers (14 mm in diameter) positioned on body landmarks based on the Motek Human Body Model marker set [42]. Marker trajectories used to capture participant motion were recorded at a sampling frequency of 100 Hz. The two multiaxial force plates (ForceLink R-Mill force plates, Motek, Netherlands), embedded within the treadmill, enabled precise readings of the ground reaction forces (GRFs) and moments at a sampling frequency of 1000 Hz. This gold standard reference system was used in the gait laboratory during the structured in-lab treadmill 2MWT.


The “silver standard” reference system consisted of body-worn IMUs manufactured by Gait Up (Physilog 6® model, Gait Up, Lausanne Switzerland). The Gait Up IMUs recorded continuous 3D accelerometer and gyroscope data at 128 Hz. Two configurations of the Gait Up system were deployed: (I) a three-sensor setup (one IMU on each foot and an additional IMU at the lower back near the center of mass) used during the structured in-lab overground and treadmill 2MWTs, and (II) a two-sensor setup, with one IMU on each foot, used during structured and unstructured real-world walking.

## Gait pipeline

A gait pipeline in biomechanics and movement analysis represents a systematic framework for evaluating human mobility [43]. It consists of an algorithmic and processing pipeline that segments raw IMU sensor signals into gait bouts, or walking segments, and detects gait events (i.e., initial and final contact), which can be used for subsequent computation of digital gait measures that measure the quantity and quality of gait. Specifically, the gait pipeline introduced here comprises the following key steps: preprocessing, alignment with gravity, task recognition, alignment with anatomical reference frame, adaptive step detection, and computation of digital gait measures assessing the quality and quantity of gait (Figure 2). However, the evaluation of the digital gait measures goes beyond the scope of the presented analysis, which focuses on the performance of the adaptive step detection algorithm in detecting initial and final contacts.

### Signal preprocessing

Initially, the data collected from the accelerometer and gyroscope were resampled to the same time grid and temporally aligned. Following this, the accelerometer's raw signal was filtered using a low-pass filter (second-order double-pass Butterworth filter with a cutoff frequency of 17 Hz [44]).

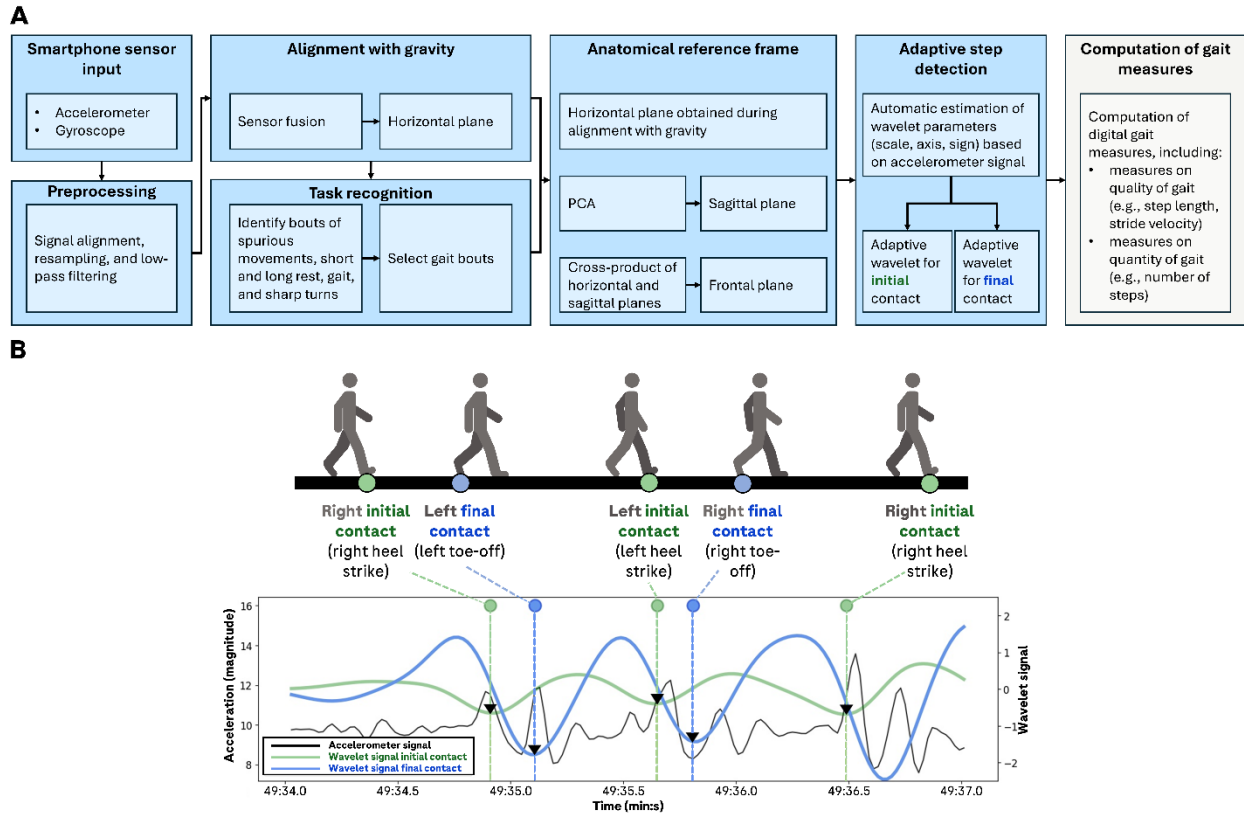

**Figure 2. Gait pipeline.** (A) The gait pipeline consists of six major steps: preprocessing, alignment with gravity, task recognition, alignment with anatomical reference frame, adaptive step detection, and gait measures calculation. The preprocessing, alignment with gravity, and alignment with anatomical reference frame steps aim to smooth and to align the accelerometer and gyroscope signals in both time and space. The task recognition step aims to identify bouts of spurious movements, short and long rests, sharp turns (i.e., turns of  $>90$  degrees), and gait bouts without sharp turns that have a duration of at least 2 seconds. Here, adaptive step detection was subsequently applied to gait bouts (the step detection algorithm was not applied to sharp turns as steps taken during sharp turns were not robustly detected by the reference measurement systems). (B) In the adaptive step detection step, parameters of the continuous wavelet transform, such as wavelet scale, wavelet axis, and wavelet sign, are automatically estimated from the accelerometer signal. Two separate wavelet transforms are applied along the axis established by the *wavelet* axis and at the scale specified by the *wavelet* scale, resulting in two wavelet signals: one represented by a green curve for detecting initial contacts (green dots) and another represented by a blue curve for detecting final contacts (blue dots). These gait events are identified in the respective wavelet transform as local minima or maxima, depending on the *wavelet* sign (here: local minima).

## Alignment with gravity

The direction of gravity in the preprocessed accelerometer and gyroscope signal, and hence the horizontal plane (Supplementary Figure S1), was determined through sensor fusion using an orientation estimation algorithm (Madgwick filter with gain parameter set to  $\beta = 0.041$  [45]). The magnetometer signal was deliberately not used in this step due to potential magnetic interference, in particular in indoor walking scenarios.

## Task recognition

Having the smartphone accelerometer and gyroscope data aligned with the direction of gravity enabled identified bouts of spurious movements, short and long rests, sharp turns, and gait (or walking). Moving and non-moving (i.e., resting) intervals were identified using both acceleration and gyroscope signals. The process involved splitting the acceleration and gyroscope signals into non-overlapping intervals. The default window length is set to 0.6 seconds, which aligns with the typical mean step duration among PwMS. However, it can be adjusted based on the specific population of interest. To identify non-moving intervals, three criteria needed to be satisfied within an interval:

- **Accelerometer magnitude:** It must fall within a defined range, specified as a “reference value  $\pm$  10% of the reference value.” The reference value is typically close to gravitational acceleration ( $g = 9.81 \text{ m/s}^2$ ).
- **Gyroscope magnitude:** This should be below a specified threshold of 0.6 rad/s, with a recommended range between 0.2 and 0.6 rad/s, where lower values more rigorously filter non-moving blocks. A 0.2 rad/s tolerance indicates static balance as a non-moving window, while 0.6 rad/s ensures adequate handling of minor movements as observed from a subset of data collected from HC, yielding satisfying results across varying smartphone locations.
- **Standard deviation of accelerometer signal:** The combined standard deviation should remain below  $0.2 \text{ m/s}^2$ . This default threshold is effective for general performance, with a recommended range from 0.05 to  $0.4 \text{ m/s}^2$  (these thresholds were determined empirically), where 0.05 identifies static balance and 0.4 allows for some movement tolerance.

Non-moving intervals less than 1 second apart were merged into a single non-moving interval. Once non-moving intervals were detected, they were classified as rests if surrounded by moving intervals, or as boundaries if they were within 2 seconds of the start or end of the recorded data. Rest intervals were further categorized based on their durations into short (less than 2 seconds) and long (equal to or longer than 2 seconds) rest intervals.

Moving intervals shorter than 2 seconds were deemed unknown, as they might include spurious movements, and were discarded in subsequent processing. The threshold of minimum gait bout duration was set at 2 seconds based on the definition of a gait bout in Kluge et al. [46], as a sequence containing at least two consecutive strides (= 4 steps). Hence, with a mean step duration is 0.5 seconds and a minimum of four steps per gait bout, the minimum gait bout duration is 2 seconds. Hence, moving intervals longer than 2 seconds were marked as possible gait bouts. After establishing the start and end points of these moving intervals (i.e., gait bouts), the segmented data underwent further analysis and step detection was applied to detect gait events.

Step detection was only applied to the gait bouts that did not contain a rest window nor sharp turns. Periods of turns were intentionally excluded as steps during turns could not be consistently enough detected with reference measurement systems. The detection of turns relied on the gyroscope signal, utilizing a well-established and commonly adopted algorithm [47]. Here, only sharp turns of  $\geq 90$  degrees were considered.

#### Verification of the moving interval as a gait segment

As gait produces a period pattern in the time domain, bouts identified as gait were verified as true gait bouts by an auto-correlation algorithm that assessed the modes of the gait bouts' auto-correlation [43].

## Alignment with anatomical reference frame

For the step detection algorithm to be agnostic to the smartphone's orientation, the smartphone sensor data must be translated to an orientation-agnostic anatomical reference frame (Supplementary Figure S1) [43]. The reference frame's horizontal plane was already obtained during alignment with gravity. The other two planes were identified as follows: The sagittal plane being parallel to the direction of movement was identified through principal component analysis (PCA) applied to the gravity-aligned accelerometer data as the direction of movement is characterized by the maximum variance (i.e., maximum movement). Finally, the frontal plan was identified by taking the cross product of the horizontal and sagittal planes. As PCA requires a sufficient amount of data—at least 3 to 4 seconds of samples—to yield reliable results [48-50], the accuracy of the identified sagittal and frontal planes was verified through a autocorrelation-based algorithm [43].

The normal vectors of these three planes were subsequently used to build the transformation matrix that defined the orientation-agnostic global coordination system. Following this transformation, the vertical axis of the smartphone sensor data was realigned with the gravitational vector and pointing upward; the antero-posterior axis was realigned with the direction of movement as determined by the PCA and pointing either forward or backward. The medio-lateral axis was then derived as the cross-product of the newly realigned vertical and antero-posterior axes, thus pointing to the left or to the right accordingly.

## Adaptive step detection

A continuous wavelet transform ('CWT'; a convolution of the acceleration data and analyzing function) identified initial contact and final contact events within the gait cycle [51]. The wavelet scale, wavelet axis (axis along which the wavelet function is applied) and wavelet sign were automatically determined based on the input acceleration and gyroscope signal [43]. This

process enhances the robustness of the pipeline by ensuring that step detection is adaptive to the characteristics of the input signals, making the step detection algorithm adaptive to different gait velocities. Once these parameters are selected, the CWT for initial and final contact detection is applied as follows:

- Integration and differentiation of wavelet axis using a Gaussian CWT, where initial contacts were identified as the timepoints of the minima or maxima based on the wavelet sign.
- The differentiated signal underwent a further CWT differentiation from which final contacts were identified as the time points of the maxima or minima based on the wavelet sign.
- The laterality of each initial contact was determined using the sign of the gyroscope signal.
- Wavelets were computed using the PyWavelets v1.8.0 package [52]

A quality check of gait events was performed to remove extraneous event detections based on physiological limits (i.e., gait events that occurred either too fast or too slow to be plausible were removed) and to verify whether the sequence of events corresponds with the gait cycle.

Furthermore, any final contacts beyond 25% of the maximum estimated stride duration were excluded. This stride duration was calculated from acceleration or gyroscope signals using an autocorrelation-based algorithm described in [43]. A 50% tolerance was applied to this estimation to obtain the maximum estimated stride duration. At the 25% mark of the gait cycle, the limb that initiated contact reaches the midstance phase, indicating that toe-off has already happened [53]. The rationale for not examining lower percentages is to ensure sufficient tolerance, thereby reducing the risk of missing any valid final contacts that may occur. Strides longer than the maximum estimated stride duration were removed.

## Statistical analysis

The analyses presented here evaluate the gait pipeline's performance in detecting gait events such as initial and final contacts. These gait events can subsequently be used to compute spatial, spatiotemporal, and temporal gait measures such as step/ stride length, stride/ step duration, stride/ step velocity as well as measures from other gait domains (e.g., measures related to gait asymmetry etc.). The evaluation of the psychometric properties of these gait measures goes, however, beyond the scope of the analyses presented in this publication and will be reported in dedicated publications.

### Gait pipeline's performance and temporal precision of gait event detection

The gait pipeline's performance was evaluated by assessing the accuracy and temporal precision with which it can identify gait events—initial and final contacts—detected by reference measurements systems (i.e., assessing the gait pipeline's technical validity). Detection performance and temporal precision were computed separately for the following groups: HC, all PwMS, PwMS with mild EDSS scores ( $\text{EDSS} \leq 3.5$ ), PwMS with moderate EDSS scores ( $\text{EDSS} 4.0 - 5.0$ ), and PwMS with severe EDSS scores ( $5.5 - 6.5$ ). Gait events detected by the smartphone were compared to those detected by the reference measurement systems (mocap and Gait Up see section “Supplementary appendix”) within a temporal tolerance window of 0.5 s, centered around the event identified by the reference tool. The 0.5-second tolerance window was selected based on its representation of a step duration and agreement with previous research [26]. Furthermore, this window also allows for some robustness in case the smartphones and reference measurement systems are not perfectly synchronized with each other, especially during structured and unstructured real-world walking, where participants were responsible for synchronizing the devices.

Gait events identified by the smartphone-based gait pipeline were classified as true positive (TP) if also identified by the reference system; false positive (FP) if not identified by the reference system; or false negative (FN) if not identified by the smartphone but identified by the reference system. If the smartphone found multiple events within the tolerance window, the event closest to the reference event was labelled a true positive. In contrast, the remaining events were marked as false positives, assuming they did not correspond to any consecutive reference event. Following the categorization of events as true positives (TP), false negatives (FN), or false positives (FP), the algorithm's performance in detecting initial and final contacts was evaluated using the following performance metrics:

$$Precision = \frac{TP}{(TP + FP)}$$

$$Recall = \frac{TP}{(TP + FN)}$$

$$F1\ score = 2 \frac{(Precision \times Recall)}{(Precision + Recall)}$$

Precision represents the proportion of true positives among all events detected by the algorithm, while recall, or sensitivity, is the proportion of true positives among all events that should have been detected. The F1 score was computed as the harmonic mean of precision and recall, providing a balanced measure of the algorithm's performance [54].

Additionally, temporal error metrics were computed that measure the temporal difference between smartphone-based and reference-based gait events (i.e., this difference is positive if the smartphone gait event occurs after the reference gait event and negative if it happens before). They thus provide insights into the accuracy and consistency of the algorithm's performance in detecting gait events. These temporal error metrics were computed across all

steps detected during a single test (i.e., all true positive gait events) and included the constant error (mean error across steps), absolute error (mean of absolute errors across steps [i.e., bias]), variable error (standard deviation of errors across steps [i.e., consistency]), total variability mean (root mean square divided by number of steps), and IQR error.

Performance and error metrics were aggregated as follows: For structured in-lab treadmill and structured in-lab overground 2MWTs—both were performed once by each participant—, performance metrics were aggregated across participants by computing the median, interquartile range (IQR), first quartile (Q1), third quartile (Q3), 5th percentile, 95th percentile, mean, and 95% confidence interval (CI) for each metric. Error metrics were aggregated by computing the mean and 95% CI across participants. In contrast, for structured and unstructured real-world walking activities, aggregation occurred in two steps as each participant completed these activities multiple times during the remote testing period. First, within-participant aggregation was performed by computing the median and IQR per participant across tests and days, resulting in one value for each metric per participant. Next, across-participant aggregation was conducted. For this, the metrics obtained from the within-participant aggregation were aggregated across participants by computing the median, IQR, Q1, Q3, 5th percentile, 95th percentile, mean, and 95% CI for performance metrics and the mean and 95% CI for error metrics. Additionally, to quantify the average of within-participants performance variability on a group level, the median of within-participants IQR was computed for the performance metrics (referred to as ‘within-subjects IQR’ [ws-IQR]).

## Effect of intrinsic, contextual and environmental gait factors on gait event detection

To further evaluate the gait pipeline’s robustness, the impact of potential confounding variables, or variable of interest, on F1 scores obtained from the structured real-world walking activity was evaluated with a Bayesian generalized linear model [55]:

$$\text{logit}(\mu[i]) = \alpha + \beta \times \overline{age} + s(\text{sex}_{idx}[i]) + d \times \sum_{j=0}^{\text{disease}_{idx}[i]-1} \delta_j + u(\text{subject}_{idx}[i]) + e(\text{environment}_{idx}[i]) + h(\text{aid}_{idx}[i])$$


Predictor variables include the global intercept ( $\alpha$ ), subject intercept to account for repeated measures within subjects ( $u$ ), age ( $\beta$ ), sex ( $s$ ), total disease effect (i.e., both disease status and disease severity;  $d$ ), contextual factors (i.e., indoors vs outdoors;  $e$ ), and walking aid usage ( $h$ ). Parameter  $\delta$ , which indicate an increment between disease status and severity, are used to constrain the effect of disease as an ordinal variable. The term “ $idx$ ” signifies index variables. Further details on the mathematical definition of the Bayesian model can be found in the section “Supplementary appendix”

## Results

### Baseline demographics and disease characteristics

A total of 35 HC and 100 PwMS were enrolled, of whom 35 HC and 96 PwMS had smartphone data available. A small subset of participants had to be excluded from the analyses due to malfunctioning reference systems, failure to use the devices during the assessments, synchronization issues between the smartphones and reference systems, incorrect identification of gait events by the reference systems, and missing or incorrect test executions. Thus, smartphone and reference measurement data were available for 32 – 35 HC and 79 – 87 PwMS across the different testing environments (Figure 3). Overall, HC and PwMS contributed towards 68 – 70 min and 158 – 172 min, respectively, of in-lab gait data, depending on the testing environment (Figure 3). A larger quantity of gait data was collected in the real-world setting as participants were instructed to complete the walking activities multiple times (i.e., once daily) rather than just once. HC contributed 548 min of gait data during the structured real-world

walking and 729 min of gait data during unstructured real-world walking. PwMS contributed towards 1422 min and 1237 min of gait data, respectively (Figure 3).

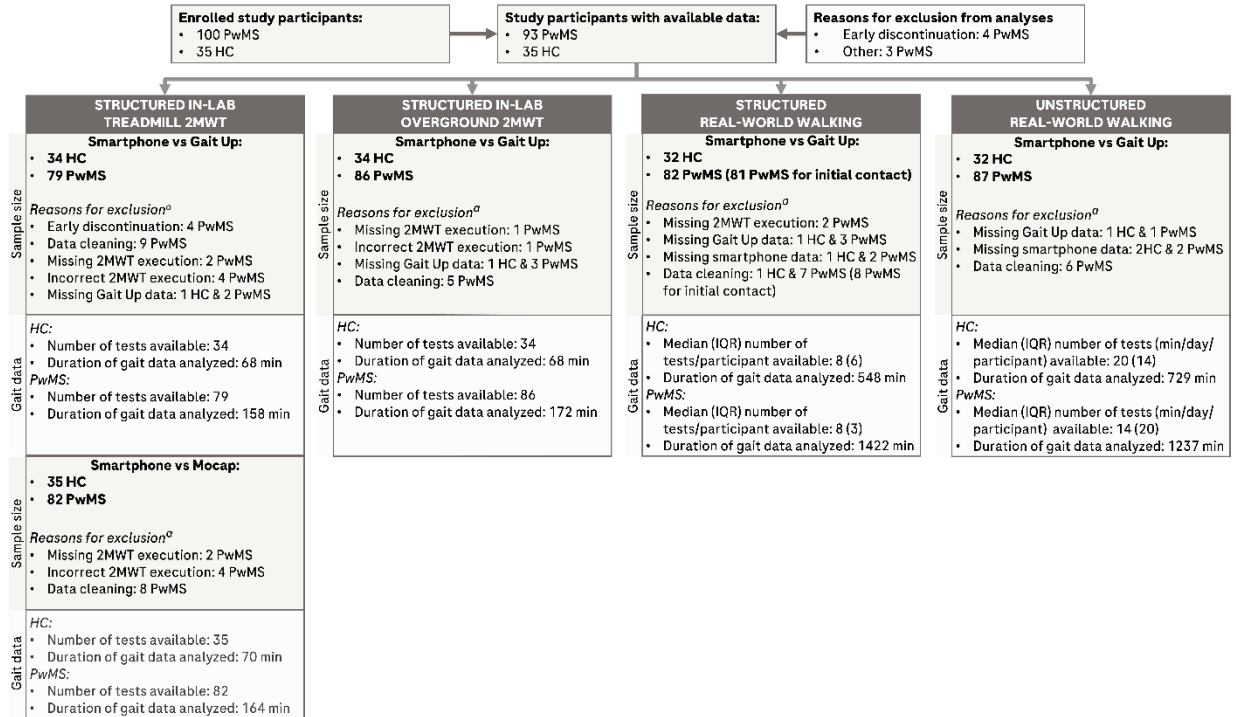

**Figure 3. Number of participants with available smartphone and reference measurement systems data and available gait data by test type.** Individual data were excluded if they did not meet the quality checks due to reference system malfunction, synchronization failure, incorrect identification of gait events by reference, and missing or incorrect test executions. <sup>a</sup>Considering all 35 HC and 100 PwMS enrolled in the study. 2MWT, Two-Minute Walk Test; HC, healthy controls; PwMS, people with multiple sclerosis.

The baseline demographics and disease characteristics are reported in Table 1. HC and PwMS were similar in age, height, weight, and body mass index. PwMS were predominantly diagnosed with relapsing-remitting disease (70% of PwMS), and the mean (SD) EDSS score among all PwMS was 4.75 (1.46). Furthermore, the PwMS cohort represented a wide range of gait impairment: T25FW times ranged between 2.9 s and 18.6 s (mean [SD]: 6.6 [3.2]), and MSWS-12 transformed total scores between 0.0 and 100.0 (mean [SD]: 45.5 [31.4]).

**Table 1. Demographics and characteristics of the study participants at baseline.**



<b>Variable</b>	<b>HC included in <math>\geq 1</math> analysis (n = 35)</b>	<b>PwMS included in <math>\geq 1</math> analysis (n = 93)</b>
Female, n (%)	19 (54)	57 (61)
Age, years, median (IQR); range	54.0 (17.5); 28.0 – 75.0	56.0 (17.0); 26.0 – 77.0
Height, m, median (IQR); range	1.7 (0.2); 1.5 – 1.9	1.7 (0.1); 1.5 – 1.9
Weight, kg, median (IQR); range	78.0 (20.0); 60.0 – 130.0	75.0 (28.0); 46.0 – 142.0
BMI, kg/cm <sup>2</sup> , median (IQR); range	26.3 (4.6); 21.5 – 34.9	25.7 (6.7); 16.9 – 45.0
MS diagnosis, n (%)		
Relapsing-remitting MS	NA	65 (70)
Secondary progressive MS	NA	15 (16)
Primary progressive MS	NA	13 (14)
Time since MS symptom onset, years, median (IQR); range	NA	18 (16.2); 3 – 53
Time since MS diagnosis, years, median (IQR); range	NA	14 (12); 1 – 39
EDSS		
Median (IQR); range	NA	5.0 (2.0); 1.5 – 6.5
0 - 3.5, n (%)	NA	23 (25)
4.0 - 5.5, n (%)	NA	35 (38)
6.0 - 6.5, n (%)	NA	35 (38)
T25FW, s, median (IQR); range	NA	5.4 (3.7); 2.9 – 18.6
MSWS-12 version 2 transformed total score, median (IQR); range	NA	40.0 (51.5); 0.0 – 100.0
MSIS-29		
Physical subscale, median (IQR); range	NA	27.5 (27.5); 0.0 - 73.8
Psychological subscale, median (IQR); range	NA	22.2 (25.0); 0.0 - 75.0

EDSS = Expanded Disability Status Scale; HC = healthy control; IQR = interquartile range; MSIS-29 = 29-item Multiple Sclerosis Impact Scale; MS = multiple sclerosis; MSWS-12 = Multiple Sclerosis Walking Scale; NA = not applicable; PwMS = people with multiple sclerosis; SD = standard deviation; T25FW = Timed 25-Foot Walk.

## Overall detection performance of gait events

### Initial contact

Structured in-lab overground and treadmill 2MWTs, performed during the second on-site visit, showed excellent median F1 scores across the six smartphone locations (HC:  $\geq 98.2\%$  for HC and  $\geq 98.5\%$  for PwMS; Figure 4a). On structured real-world walking, median (IQR) F1 scores consistently remained high, with HC at 100% (0.3%) and PwMS at 99.5% (1.9%). The lowest median (IQR) F1 score was observed in the severe cohort, which recorded 98.5% (6.0%) (Figure 4a). In contrast, unstructured real-world walking activities involving free walking yielded a median (IQR) F1 score of 97.8% (2.6%) in HC and 94.4% (6.2%) in the MS cohort with a trend of higher F1 scores for less impaired patients (Figure 4a). Median (IQR) F1 scores varied from 91.8% (8.1%) for severe PwMS to 96.0% (3.6%) for mild PwMS. Thus, although there was a small performance drop in unstructured real-world walking compared to other test types (i.e., testing environments), the accuracy of the step detection remained very high.

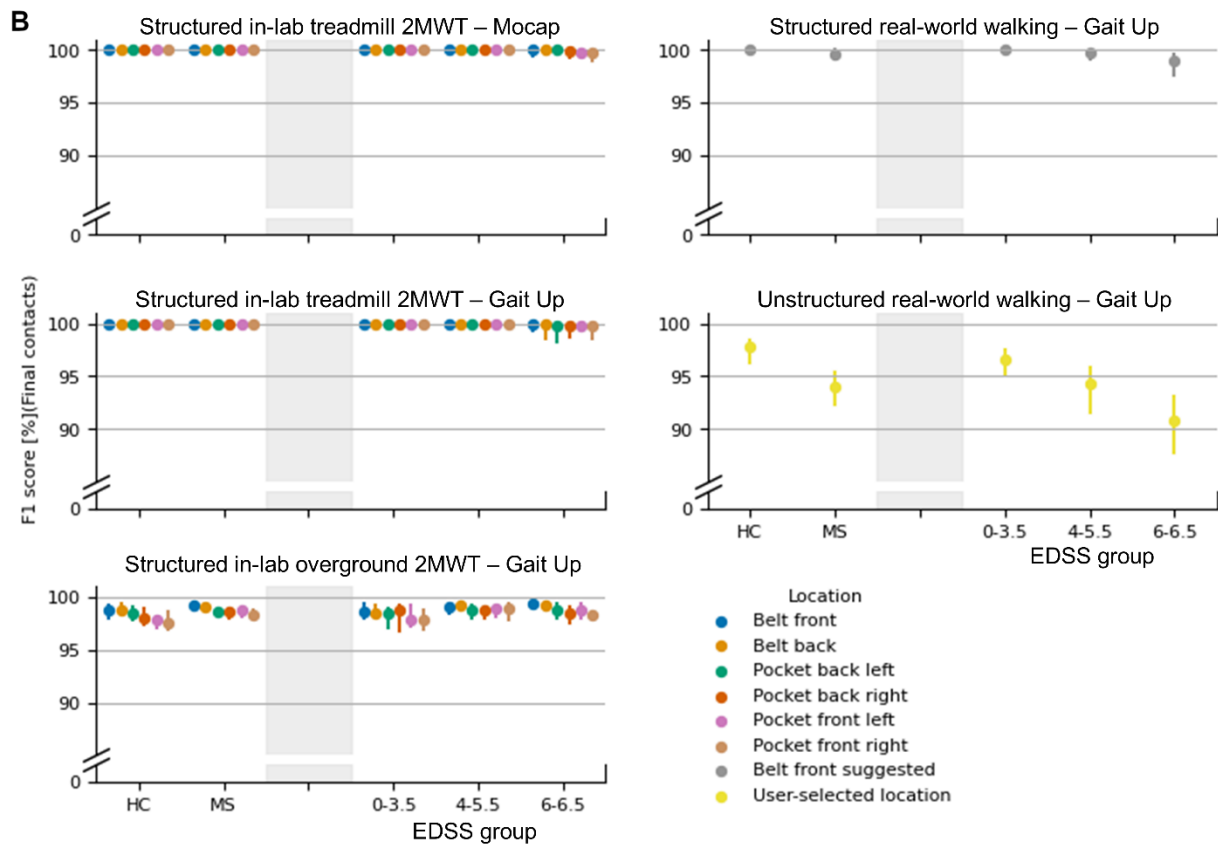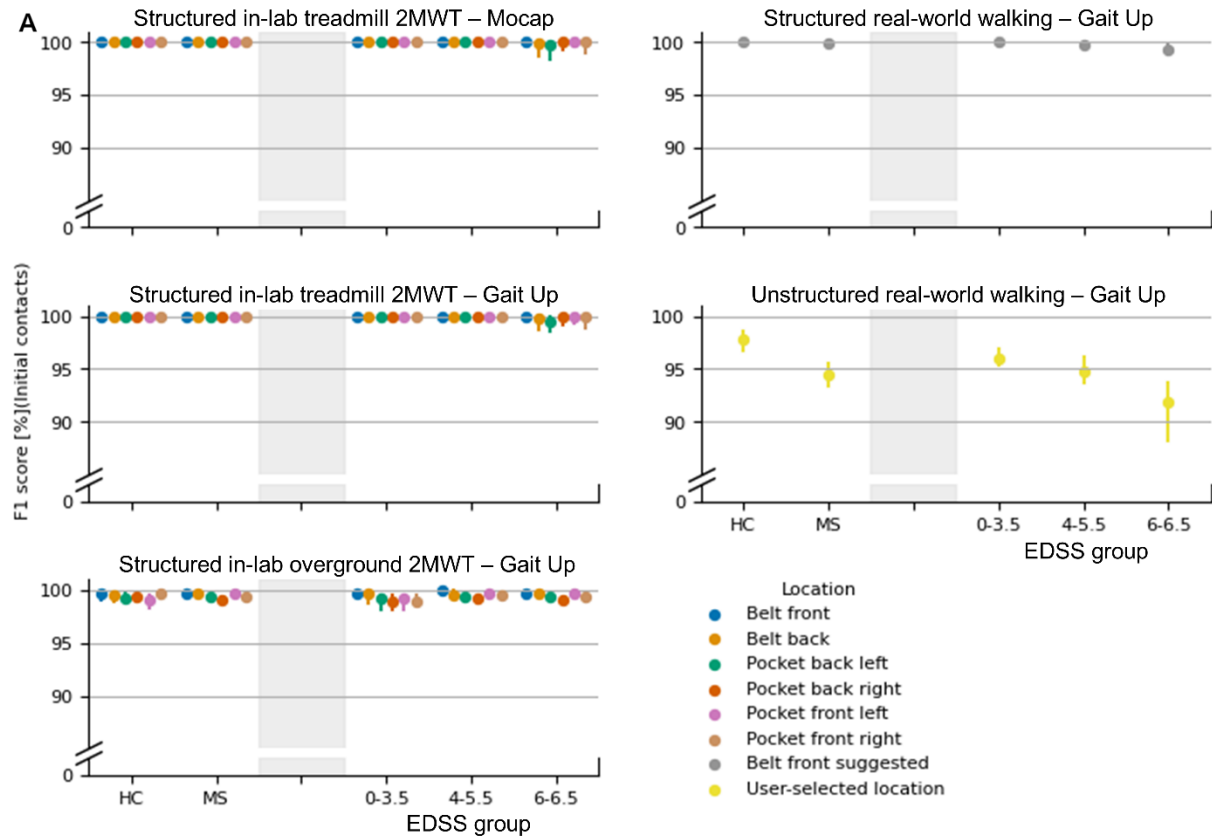

**Figure 4. Median F1 score (IQR) by smartphone wear location, cohort, and test setting.** (A) F1 scores are shown for initial contacts. (B) F1 scores are shown for final contacts. In both panels, note the break in the y-axes (F1 scores) between 5% and 85%, which was introduced for better readability.

Precision and recall scores were also high (Supplementary Table S1). Across test types and smartphone wear locations, median precision scores were  $\geq 97.4\%$  in HC and  $\geq 93.0\%$  in PwMS, whereas median recall scores were  $\geq 98.7\%$  in HC and  $\geq 96.4\%$  in PwMS.

Detailed statistics—including medians, means, quantiles, and interquartile ranges for the F1, precision, and recall scores—categorized by smartphone location, cohort, and setting are reported in Supplementary Table S1.

## Final contact

Similar findings were obtained for the final contacts. Both structured in-lab overground and treadmill 2MWTs showed very high median F1 scores across smartphone wear locations (HC:  $\geq 96.5\%$ ; PwMS:  $\geq 97.7\%$ ; Figure 4b), indicating a 2% drop in median F1 scores on the structured in-lab overground 2MWT compared with initial contacts (final vs initial contact:  $\geq 96.5\%$  vs  $\geq 98.2\%$  for HC;  $\geq 97.7\%$  vs  $\geq 98.5\%$  for PwMS). Consistently high median (IQR) F1 scores were also observed on structured real-world walking (100% [0.2%] for HC; 99.4% [2.5%] for PwMS) (Figure 4b). Across disease severity groups, F1 scores varied from 98.1% (6.5%) in severe PwMS to 100% (0.1%) in mild PwMS. High median (IQR) F1 scores were also obtained on unstructured real-world walking (HC: 97.8% [2.8%]; PwMS: 94.0% [6.5%] PwMS; Figure 4b). A trend of higher F1 in less impaired PwMS was observed, with F1 scores varying between 90.9% (9.1%) in severe PwMS to 96.5% (4.1%) in mild PwMS.

Precision and recall scores were equally high across test types and smartphone wear locations. Median precision scores were  $\geq 94.9$  for HC and  $\geq 93.6\%$  for PwMS, whereas median recall scores ranged were  $\geq 98.1\%$  for HC and 95.4% for PwMS.

Detailed statistics—including medians, means, quantiles, and interquartile ranges for the F1 score, precision score, and recall score—categorized by smartphone location, cohort, and testing environment are reported in Supplementary Table S2.

## Temporal precision of gait event detection

### Initial contact


The median errors recorded during the remote setting were comparable to those observed in-lab settings and between HC and PwMS (Figure 5a). Both structured in-lab overground and treadmill 2MWTs indicated that the median errors across the six smartphone locations in HC and PwMS ranged between 0.02 s to 0.08 s (Figure 5a). Notably, smartphone locations on the front belt, back belt, and front pockets consistently yielded lower median errors, ranging from 0.02 s to 0.07 s, compared to those in back pockets, which exhibited median errors between 0.05 s and 0.08 s. Few smartphone wear locations such as the front belt, back belt and front pockets showed a mild trend of larger temporal errors being associated with higher EDSS scores; however, no such trend was observed for the back pockets (Figure 5a)

The median (IQR) errors for structured and unstructured real-world walking activities were similar, although unstructured real-world walking tended to produce slightly larger errors (Figure 5a). On structured real-world walking, median (IQR) errors were 0.05 s (0.02 s) in HC and 0.06 s (0.03 s) in PwMS. In contrast, for unstructured real-world walking activities, median (IQR) errors were 0.06 s (0.02 s) in HC and 0.07 s (0.02 s) in PwMS. The median errors remained stable with

increasing EDSS scores during both structured and unstructured real-world walking (Figure 5a). Detailed statistics on temporal errors for initial contact are reported in Supplementary Table S3.

## Final contact

Median errors varied across smartphone locations in the in-lab setting, with smartphone placements on the front belt, back belt, and front pockets consistently demonstrating lower median errors in HC, but not in PwMS, (median temporal error range: 0.02 s – 0.06 s in HC; 0.01 s to s 0.05 in PwMS) when compared to those situated in back pockets (median temporal error range: 0.04 s – 0.08 s in HC; 0.01 s – 0.05 s in PwMS) (Fig 5b). The median errors for structured and unstructured real-world walking activities were similar; however, the unstructured walking demonstrated greater variability in these errors. On structured real-world walking, median (IQR) errors were 0.08 s (0.03 s) in HC and 0.05 s (0.07 s) in PwMS. On unstructured real-world walking activities, median (IQR) errors were 0.07 s (0.02 s) in HC and 0.06 s (0.06 s) in PwMS. Furthermore, the median errors were slightly lower in PwMS than in HC across varying environmental settings. As the severity of the disease progressed, the median errors tended to improve in both remote and in-lab settings. However, the variability in these errors also increased with higher EDSS scores, indicating that more severe cases may lead to less consistent results. Detailed statistics on temporal errors for final contact are reported in Supplementary Table S4.

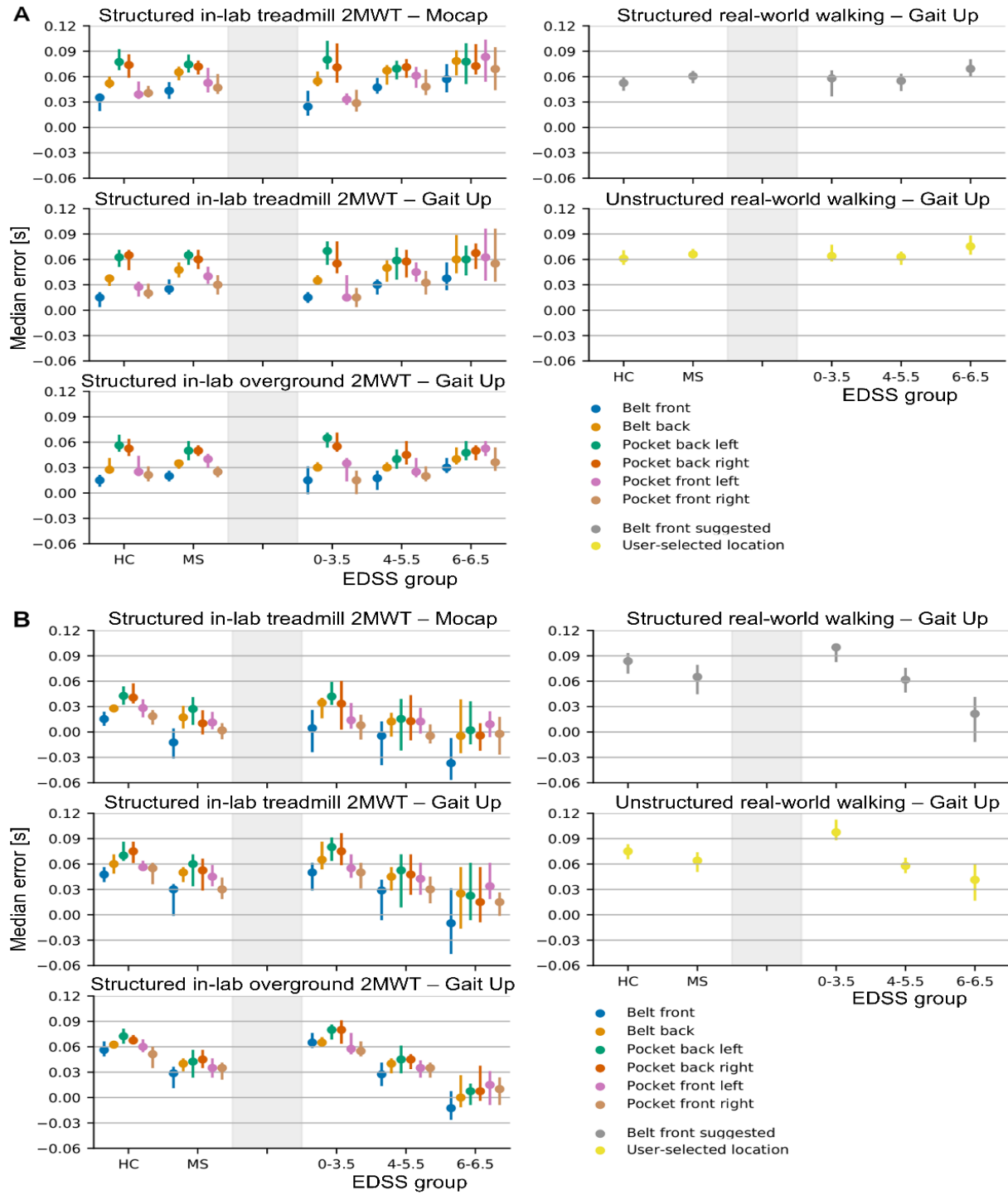

**Figure 5. Median temporal error (mean [95% CI]) by smartphone wear location, cohort, and test type.** Median temporal errors for initial contacts (A) and final contacts (B) were aggregated across participants by computing the mean and 95% CI across participants.


## Effect of intrinsic, contextual and environmental factors on gait event detection

### Initial contact

The impact of the examined factors (i.e., age and gender, use of walking aid, indoors versus outdoors setting, and disease status and severity) on F1 score is minimal. None of the factors reach statistical significance, as evidenced by distributions crossing zero (Figure 6a). While statistically non-significant, the strongest differences were identified for disease status and severity and the comparison between the indoors and outdoors settings. In both scenarios, a decrease of about 1.7-2% in the median F1 score was observed for disease severity, as well as when participants walked indoors compared to their outdoor walking. More detailed statistics, including parameter distributions, are presented in Supplementary Table S5.

### Final contact

Results for the final contact are highly consistent with those reported for the initial contacts and are summarized in Figure 6b and Supplementary Table S6.

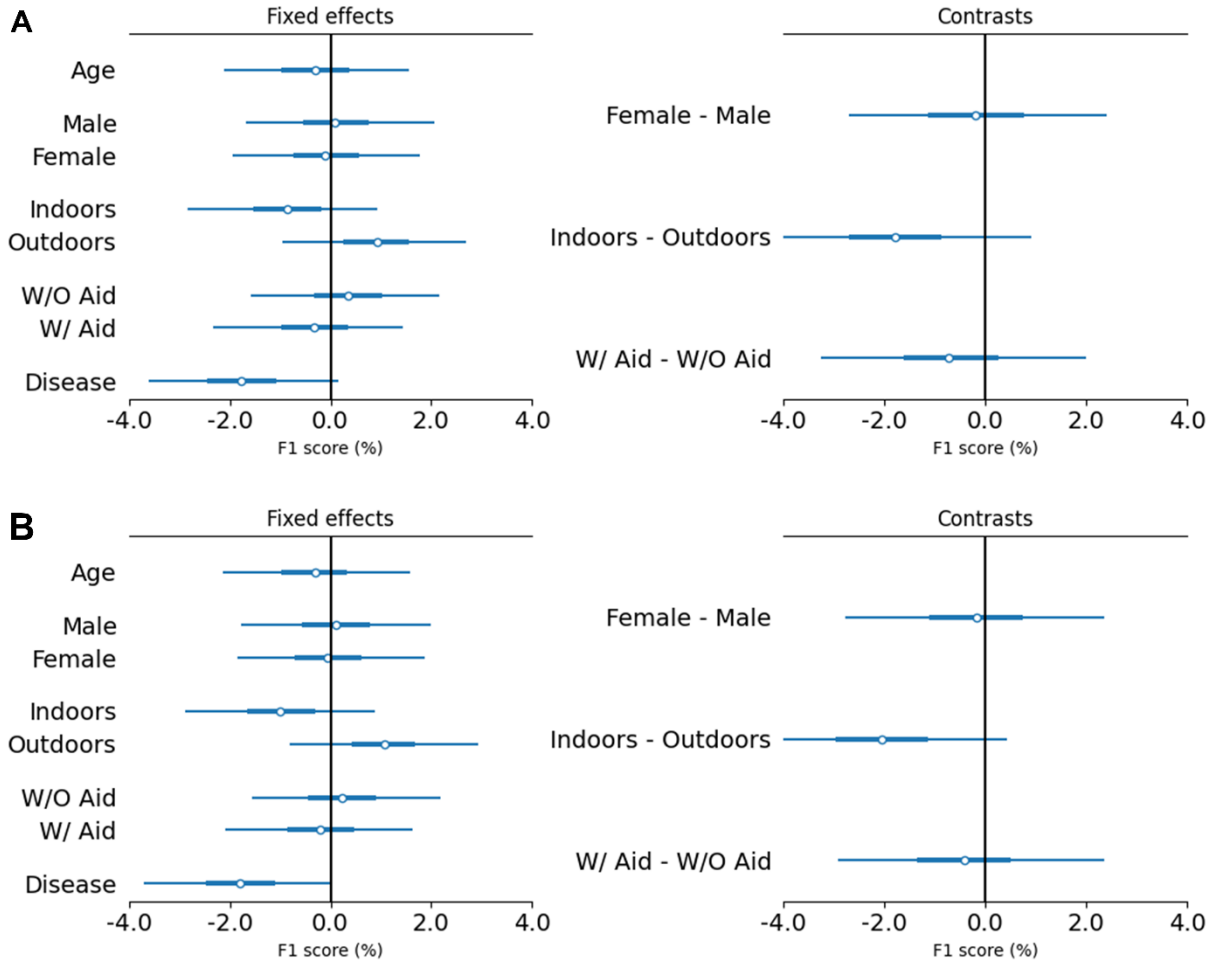

**Figure 6. Parameter estimates of Bayesian model examining the influence of different factors on the F1 score.** The forest plots depict the credible intervals of both the fixed effects and contrast of the Bayesian model for initial contacts (A) and final contacts (B). The central points are the estimated posterior median, the thick lines are the IQR range (0.25-0.75), and the thin lines represent the highest density intervals (94%). The forest plots depict both the fixed effects and contrast of the Bayesian model for initial contacts (A) and final contacts (B). Contrast depicts the difference in the fixed effects between related factors such as sex (female vs male), test setting (indoors vs outdoors), are use of walking aids (with vs without walking aids). W/ Aid, with walking aids or orthotics; W/O Aid, without walking aids and/ or orthotics.

## Discussion

This work introduces a gait pipeline with adaptive step detection that accurately identifies gait events across different walking activities in both HC and PwMS with varying disease severity. The results demonstrated the gait pipeline's ability to detect both initial and final contacts with high precision and high recall across multiple environments (controlled gait laboratory setting vs real-world scenarios), different smartphone sensor positions, and different conditions (e.g.,

indoor vs outdoor, with or without use of walking aids). These results confirm the robust technical validity of the pipeline and its potential for real-world unsupervised gait assessment. Importantly, this gait pipeline addresses one of the key challenges in advancing real-world gait analysis, the translation of wearable algorithms designed for and validated in controlled indoor settings to real-world unsupervised settings [36, 56].

## Overall performance of gait events detection

The gait pipeline demonstrated a precision in identifying initial contacts in unstructured real-world walking that is comparable to the best-performing algorithms for detecting initial contacts [33, 51, 57-60], recently evaluated by Micó-Amigo et al. [61] (mean [95% CI] precision: 91% [90%, 93%] vs 92% [91%, 93%]). Additionally, the gait pipeline achieved higher recall than the best-performing algorithm evaluated by Micó-Amigo et al. [61] (mean [95%CI] recall: 94% [93%, 96%] vs 82% [81% - 84%]). This greater performance compared with the work by Micó-Amigo et al.[61] was achieved even though the technical validation was conducted in more challenging conditions that could have potentially contributed to lower precision and recall metrics: (I) longer data collection periods (2 weeks vs 2.5 hours) and resulting larger number of gait bouts, (II) a larger cohort of more severely impaired PwMS who are more likely to use walking aids or orthotics (mean [SD] baseline EDSS: 4.75 [1.46] vs 3.5 [1.7]), (III) additional sensor positions (6 devices vs 1 device), and (IV) less controlled settings (both in-lab and real-world structured and unstructured walking activities vs specific tasks in real-world setting). Consequently, the data used for the technical validation of the gait pipeline captured a wider range of real-world walking variability, underscoring its robustness.

Previous research has highlighted the challenge of comparing results obtained in controlled settings with those from unsupervised, real-world settings due to differences in patient behaviors

and movement complexity in natural surroundings [36, 62]. Remarkably, our gait pipeline consistently demonstrated high performance in both in-lab and real-world settings, regardless of the type of population (pathological and non-pathological) and disease stage. As can be expected, F1 scores were slightly lower in the unstructured real-world walking activities due to a greater number of shorter gait bouts lasting less than 10s as well as a greater proportion of transitional phases of gait (e.g., initiation and termination of gait), turns, and changes in walking speed, all of which may negatively impact step detection performance. However, the reduction in F1 scores was minimal (approximately 2% points for HC and 5% points for PwMS), and the F1 scores remained higher than those reported for previously published algorithms. The pipeline's improved performance can be explained with the novel adaptive step detection algorithm, which uses an adaptive wavelet to enhance the detection of initial and final contacts in different settings, and the robust anatomical reference system identification. The wavelet scale, sign, and the axis on which it is applied are determined based on the input acceleration signal. This adaptability allows the pipeline to increase its robustness in tackling several challenges, including significant variability in signal characteristics, intrinsic gait factors, such as population differences (e.g., healthy adults versus patients with various conditions and stages of disease) and variations in walking speed (e.g., slow versus fast), contextual factors like the location of walks (e.g., indoors vs. outdoors) and the use of walking aids, and environmental factors, which encompass differences between laboratory, clinical, and real-world settings. Together, these results suggest the potential utility of the gait pipeline for future assessments of pathological gait in natural and unobtrusive settings.

## Temporal precision of gait event detection

Temporal errors identifying initial contact were consistent with those reported in the existing literature. These errors remained relatively stable across various settings, including in-lab and real-world settings, structured and unstructured walking activities, different smartphone wear

locations, and varying levels of disability. Differences in temporal errors were noted in specific cases. For instance, smartphones placed in the back pocket tended to exhibit slightly larger temporal errors compared to other smartphone wear locations, likely due to increased inference from the soft tissue artifacts. Furthermore, real-world walking activities tended to show larger temporal errors than in-lab walking activities. Nonetheless, temporal errors for real-world walking fell within the same range as previously reported for CWT-based algorithms applied to in-lab gait data [33, 61, 63].

All these differences in temporal errors lie within the margin of uncertainty that should be considered when interpreting time differences in gait events. This uncertainty stems from the measurement error inherent in reference measurement systems and can be as large as 0.04 s. Specifically, the event detection algorithms incorporated in the Gait Offline Analysis Tool (GOAT) software package—employed to establish the reference initial and final contacts referenced in this paper—(see section “Supplementary appendix” for further details) reported temporal errors as large as 0.03 s in impaired population such PwMS or people with stroke [64, 65].

Most previous MS studies reported temporal errors only for initial contact and typically restricted to a single sensor location, the lower back. Assessments of temporal errors in identifying final contact have been mostly conducted in other patient populations, such as people with Parkinson’s disease [33]. But as gait is differently impacted in PwMS compared with people with Parkinson’s disease, a direct comparison of temporal errors across these two populations is difficult. Here, we expanded on previous MS studies by evaluating multiple smartphone wear locations and, crucially, also by assessing temporal errors in identifying final contact. We found that temporal errors for final contact were in the same range of magnitude as those for initial contact. Significantly, the constant positive bias we observed is not expected to have any impact on the computation of such gait measures as cadence, step symmetry, and gait variability. This

suggests the gait pipeline detects both types of gait events with a high degree of temporal accuracy.

## Effect of intrinsic, contextual and environmental factors on gait event detection

While minor variations were observed, the overall detection performance remained robust and reliable across different factors, including disease severity—the pipeline showed good performance in detecting steps in PwMS with varying levels of MS-related disability (mild, moderate, and severe)—and environmental conditions—the pipeline maintained high performance in both indoor and outdoor settings.

After examining the information from the gait diary on how the structured real-world walking was performed, a Bayesian analysis found that factors like indoor vs. outdoor setting, use of walking aids, age, gender, population type, and disease stage had minimal and statistically non-significant impact on the F1 score. Age and gender did not affect the percentage difference in F1 score for initial and final contacts. The median F1 score was about 2% lower in indoor compared to outdoor settings for initial and final contacts. This small difference may be due to the more confined nature of indoor settings [30], potentially leading to a higher frequency of shorter gait bouts and an increased presence of transitional phases of gait and turns. Using walking aids can also cause changes in walking patterns, which in turn may directly impact the pipeline's performance [66, 67]. The performance of our gait pipeline was marginally reduced by 0.7% for initial contacts and 0.4% for final contacts in PwMS who used walking aids. Nonetheless, this reduction in performance did not reach statistical significance, indicating that while there is a notable impact, the use of walking aids does not substantially compromise the overall performance of the gait pipeline. Moreover, the pipeline maintains high performance across

varying levels of MS-related disability, or MS disease severity. A slight drop of approximately 2% in the median F1 score as disability increased; however, this change did not reach statistical significance, reaffirming the pipeline's robustness regardless of severity.

These findings indicate that the gait pipeline's gait performance remained consistent regardless of intrinsic gait factors as well as contextual and environmental influences. Even with a slight drop in performance, the pipeline performed to a very high level of accuracy. The demographics and anthropometric characteristics of the enrolled PwMS exhibit considerable diversity. Specifically, PwMS with various clinical phenotypes of MS (i.e., relapsing-remitting MS, primary progressive MS, and secondary progressive MS), along with varying levels of disability (i.e., mild, moderate, and severe cohorts), were included in the analyses. This heterogeneity within the study population establishes a robust framework for objectively assessing the performance of the gait pipeline in MS. Additionally, it enhances the potential for generalizing the findings to a broader MS population. Furthermore, the pipeline has the added benefit of not requiring any knowledge of the exact sensor location and orientation of the smartphones, contrary to many previously validated algorithms [32, 33, 51, 61, 62]. Consequently, this gait pipeline has less stringent requirements for future study participants or patients on how to carry the smartphones (or other IMU sensors) while performing the gait assessments.

## Limitations

This study comes with several limitations. First, the gait pipeline's performance in detecting steps during turns could not be evaluated due to the absence of motion capture capabilities during turns—the motion capture system was only available for the structured in-lab treadmill 2MWT which did not involve turns—and due to challenges in consistently detecting steps during turns with other reference measurement systems. Second, there were challenges related to

participant adherence, data loss, and synchronization problems. Nonetheless, the sample size available for the analyses was large enough (79 – 87 of the 100 PwMS enrolled and 32 – 35 of the 35 HC enrolled, depending on the analysis) to support the robustness of the gait pipeline. Third, the real-world walking activities were performed with a single smartphone device. Consequently, a within-subject comparison of different smartphone wear locations was not possible for these activities.

## Conclusions

This study demonstrates the technical validity and robustness of a gait pipeline that uses an adaptive step detection algorithm for use in PwMS, highlighting its potential for real-world gait assessment. Crucially, the pipeline performed with high accuracy across different test settings, varying levels of disease severities, different demographic characteristics, and different environmental conditions. These findings expand on previous research by not only validating the detection of initial contacts but also the detection of final contacts. Such accurate step detection is essential for reliably estimating digital gait measures such as cadence, step symmetry, and gait variability. It is, therefore, expected that the findings will advance our understanding of gait analysis and contribute to this field's existing knowledge. By developing a new and robust gait pipeline that offers flexibility to the setting in which gait data are captured and where the smartphone is worn could pave the way for an assessment of pathological gait that is not only accurate and reliable but also easily deployable in the real world.

# References

- [1] A. J. Thompson, S. E. Baranzini, J. Geurts, B. Hemmer, and O. Ciccarelli, "Multiple sclerosis," *The Lancet*, vol. 391, no. 10130, pp. 1622-1636, 2018.
- [2] L. Comber, R. Galvin, and S. Coote, "Gait deficits in people with multiple sclerosis: A systematic review and meta-analysis," *Gait Posture*, vol. 51, pp. 25-35, Jan, 2017.
- [3] N. G. LaRocca, "Impact of Walking Impairment in Multiple Sclerosis," *The Patient: Patient-Centered Outcomes Research*, vol. 4, no. 3, pp. 189-201, 2011/09/01, 2011.
- [4] M. Yildiz, "The impact of slower walking speed on activities of daily living in patients with multiple sclerosis," *International Journal of Clinical Practice*, vol. 66, no. 11, pp. 1088-1094, 2012/11/01, 2012.
- [5] C. Heesen, J. Böhm, C. Reich, J. Kasper, M. Goebel, and S. M. Gold, "Patient perception of bodily functions in multiple sclerosis: gait and visual function are the most valuable," *Mult Scler*, vol. 14, no. 7, pp. 988-91, Aug, 2008.
- [6] G. Kobelt, A. Thompson, J. Berg, M. Gannedahl, and J. Eriksson, "New insights into the burden and costs of multiple sclerosis in Europe," *Mult Scler*, vol. 23, no. 8, pp. 1123-1136, Jul, 2017.
- [7] D. A. Scalzitti, K. J. Harwood, J. R. Maring, S. J. Leach, E. A. Ruckert, and E. Costello, "Validation of the 2-Minute Walk Test with the 6-Minute Walk Test and Other Functional Measures in Persons with Multiple Sclerosis," *Int J MS Care*, vol. 20, no. 4, pp. 158-163, Jul-Aug, 2018.
- [8] R. W. Motl, J. A. Cohen, R. Benedict, G. Phillips, N. LaRocca, L. D. Hudson, and R. Rudick, "Validity of the timed 25-foot walk as an ambulatory performance outcome measure for multiple sclerosis," *Mult Scler*, vol. 23, no. 5, pp. 704-710, Apr, 2017.
- [9] E. Sebastião, B. M. Sandroff, Y. C. Learmonth, and R. W. Motl, "Validity of the Timed Up and Go Test as a Measure of Functional Mobility in Persons With Multiple Sclerosis," *Arch Phys Med Rehabil*, vol. 97, no. 7, pp. 1072-7, Jul, 2016.
- [10] L. Kappos. "Neurostatus scoring definitions: Version 04/10.2. Available via: [https://www.neurostatus.net/media/specimen/Definitions\\_0410-2\\_s.pdf](https://www.neurostatus.net/media/specimen/Definitions_0410-2_s.pdf) [Accessed January 2025]."
- [11] A. G. Skjerbæk, L. G. Hvid, F. Boesen, L. Taul-Madsen, E. Stenager, and U. Dalgas, "Psychometric measurement properties and reference values of the six-spot step test, the six-minute walk test, the 25-foot walk test, and the 12-item multiple sclerosis walking scale in people with multiple sclerosis," *Multiple Sclerosis and Related Disorders*, vol. 94, pp. 106242, 2025/02/01/, 2025.
- [12] M. E. L. van den Berg, C. J. Barr, J. V. McLoughlin, and M. Crotty, "Effect of walking on sand on gait kinematics in individuals with multiple sclerosis," *Multiple Sclerosis and Related Disorders*, vol. 16, pp. 15-21, 2017/08/01/, 2017.
- [13] L. Angelini, W. Hodgkinson, C. Smith, J. M. Dodd, B. Sharrack, C. Mazzà, and D. Paling, "Wearable sensors can reliably quantify gait alterations associated with disability in people with progressive multiple sclerosis in a clinical setting," *J Neurol*, vol. 267, no. 10, pp. 2897-2909, Oct, 2020.
- [14] S. Shema-Shiratzky, E. Gazit, R. Sun, K. Regev, A. Karni, J. J. Sosnoff, T. Herman, A. Mirelman, and J. M. Hausdorff, "Deterioration of specific aspects of gait during the instrumented 6-min walk test among people with multiple sclerosis," *J Neurol*, vol. 266, no. 12, pp. 3022-3030, Dec, 2019.
- [15] S. Alexander, G. Peryer, E. Gray, F. Barkhof, and J. Chataway, "Wearable technologies to measure clinical outcomes in multiple sclerosis: A scoping review," *Mult Scler*, vol. 27, no. 11, pp. 1643-1656, Oct, 2021.

- [16] L. Angelini, I. Carpinella, D. Cattaneo, M. Ferrarin, E. Gervasoni, B. Sharrack, D. Paling, K. P. S. Nair, and C. Mazzà, "Is a Wearable Sensor-Based Characterisation of Gait Robust Enough to Overcome Differences Between Measurement Protocols? A Multi-Centric Pragmatic Study in Patients with Multiple Sclerosis," *Sensors (Basel)*, vol. 20, no. 1, Dec 21, 2019.
- [17] R. I. Spain, R. J. St George, A. Salarian, M. Mancini, J. M. Wagner, F. B. Horak, and D. Bourdette, "Body-worn motion sensors detect balance and gait deficits in people with multiple sclerosis who have normal walking speed," *Gait Posture*, vol. 35, no. 4, pp. 573-8, Apr, 2012.
- [18] J. J. Craig, A. P. Bruetsch, S. G. Lynch, F. B. Horak, and J. M. Huisinga, "Instrumented balance and walking assessments in persons with multiple sclerosis show strong test-retest reliability," *J Neuroeng Rehabil*, vol. 14, no. 1, pp. 43, May 22, 2017.
- [19] M. De Angelis, L. Lavorgna, A. Carotenuto, M. Petruzzo, R. Lanzillo, V. Brescia Morra, and M. Moccia, "Digital Technology in Clinical Trials for Multiple Sclerosis: Systematic Review," *J Clin Med*, vol. 10, no. 11, May 26, 2021.
- [20] E. Walter, M. Traunfellner, F. Meyer, C. Enzinger, M. Guger, C. Bsteh, P. Altmann, H. Hegen, C. Goger, and V. Mikl, "Cost-effectiveness of the Floodlight(®) MS app in Austria. Unlocking the mystery of costs and outcomes of a digital health application for patients with multiple sclerosis," *Digit Health*, vol. 11, pp. 20552076251314550, Jan-Dec, 2025.
- [21] L. Midaglia, P. Mulero, X. Montalban, J. Graves, S. L. Hauser, L. Julian, M. Baker, J. Schadrack, C. Gossens, A. Scotland, F. Lipsmeier, J. van Beek, C. Bernasconi, S. Belachew, and M. Lindemann, "Adherence and Satisfaction of Smartphone- and Smartwatch-Based Remote Active Testing and Passive Monitoring in People With Multiple Sclerosis: Nonrandomized Interventional Feasibility Study," *J Med Internet Res*, vol. 21, no. 8, pp. e14863, 2019/08/30, 2019.
- [22] A. P. Creagh, C. Simillion, A. K. Bourke, A. Scotland, F. Lipsmeier, C. Bernasconi, J. v. Beek, M. Baker, C. Gossens, M. Lindemann, and M. D. Vos, "Smartphone- and Smartwatch-Based Remote Characterisation of Ambulation in Multiple Sclerosis During the Two-Minute Walk Test," *IEEE Journal of Biomedical and Health Informatics*, vol. 25, no. 3, pp. 838-849, 2021.
- [23] X. Montalban, J. Graves, L. Midaglia, P. Mulero, L. Julian, M. Baker, J. Schadrack, C. Gossens, M. Ganzetti, A. Scotland, F. Lipsmeier, J. van Beek, C. Bernasconi, S. Belachew, M. Lindemann, and S. L. Hauser, "A smartphone sensor-based digital outcome assessment of multiple sclerosis," *Multiple Sclerosis Journal*, vol. 28, no. 4, pp. 654-664, 2022.
- [24] P. van Oirschot, M. Heerings, K. Wendrich, B. den Teuling, F. Dorssers, R. van Ee, M. B. Martens, and P. J. Jongen, "A Two-Minute Walking Test With a Smartphone App for Persons With Multiple Sclerosis: Validation Study," *JMIR Form Res*, vol. 5, no. 11, pp. e29128, Nov 17, 2021.
- [25] M. Scaramozza, A. Ruet, P. A. Chiesa, L. Ahamada, E. Bartholomé, L. Carment, J. Charre-Morin, G. Cosne, L. Diouf, C. C. Guo, A. Juraver, C. M. Kanzler, A. Karatsidis, C. Mazzà, J. Penalver-Andres, M. Ruiz, A. Saubusse, G. Simoneau, A. Scotland, Z. Sun, M. Tang, J. van Beek, L. Zajac, S. Belachew, B. Brochet, and N. Campbell, "Sensor-Derived Measures of Motor and Cognitive Functions in People With Multiple Sclerosis Using Unsupervised Smartphone-Based Assessments: Proof-of-Concept Study," *JMIR Form Res*, vol. 8, pp. e60673, Nov 8, 2024.
- [26] R. Romijnders, F. Salis, C. Hansen, A. Küderle, A. Paraschiv-Ionescu, A. Cereatti, L. Alcock, K. Aminian, C. Becker, S. Bertuletti, T. Bonci, P. Brown, E. Buckley, A. Cantu, A. E. Carsin, M. Caruso, B. Caulfield, L. Chiari, I. D'Ascanio, S. Del Din, B. Eskofier, S. J. Fernstad, M. S. Fröhlich, J. Garcia Aymerich, E. Gazit, J. M. Hausdorff, H. Hiden, E. Hume, A. Keogh, C. Kirk, F. Kluge, S. Koch, C. Mazzà, D. Megaritis, E. Micó-Amigo, A. Müller, L. Palmerini, L. Rochester, L. Schwickert, K. Scott, B. Sharrack, D. Singleton, A.

- Soltani, M. Ullrich, B. Vereijken, I. Vogiatzis, A. Yarnall, G. Schmidt, and W. Maetzler, "Ecological validity of a deep learning algorithm to detect gait events from real-life walking bouts in mobility-limiting diseases," *Front Neurol*, vol. 14, pp. 1247532, 2023.
- [27] F. A. Storm, C. J. Buckley, and C. Mazzà, "Gait event detection in laboratory and real life settings: Accuracy of ankle and waist sensor based methods," *Gait & Posture*, vol. 50, pp. 42-46, 2016/10/01/, 2016.
- [28] H.-h. Lee, S. Choi, and M.-j. Lee, "Step Detection Robust against the Dynamics of Smartphones," *Sensors*, vol. 15, no. 10, pp. 27230-27250, 2015.
- [29] M. S. Pan, and H. W. Lin, "A Step Counting Algorithm for Smartphone Users: Design and Implementation," *IEEE Sensors Journal*, vol. 15, no. 4, pp. 2296-2305, 2015.
- [30] F. A. Storm, B. W. Heller, and C. Mazzà, "Step detection and activity recognition accuracy of seven physical activity monitors," *PLoS One*, vol. 10, no. 3, pp. e0118723, 2015.
- [31] M. Tietsch, A. Muaremi, I. Clay, F. Kluge, H. Hoefling, M. Ullrich, A. Küderle, B. M. Eskofier, and A. Müller, "Robust Step Detection from Different Waist-Worn Sensor Positions: Implications for Clinical Studies," *Digit Biomark*, vol. 4, no. Suppl 1, pp. 50-58, Winter, 2020.
- [32] R. Romijnders, E. Warmerdam, C. Hansen, J. Welzel, G. Schmidt, and W. Maetzler, "Validation of IMU-based gait event detection during curved walking and turning in older adults and Parkinson's Disease patients," *J Neuroeng Rehabil*, vol. 18, no. 1, pp. 28, Feb 6, 2021.
- [33] M. H. Pham, M. Elshehabi, L. Haertner, S. Del Din, K. Srulijes, T. Heger, M. Synofzik, M. A. Hobert, G. S. Faber, C. Hansen, D. Salkovic, J. J. Ferreira, D. Berg, Á. Sanchez-Ferro, J. H. van Dieën, C. Becker, L. Rochester, G. Schmidt, and W. Maetzler, "Validation of a Step Detection Algorithm during Straight Walking and Turning in Patients with Parkinson's Disease and Older Adults Using an Inertial Measurement Unit at the Lower Back," *Front Neurol*, vol. 8, pp. 457, 2017.
- [34] M. D. Czech, D. Psaltos, H. Zhang, T. Adamusiak, M. Calicchio, A. Kelekar, A. Messere, K. R. A. Van Dijk, V. Ramos, C. Demanuele, X. Cai, M. Santamaria, S. Patel, and F. I. Karahanoglu, "Age and environment-related differences in gait in healthy adults using wearables," *NPJ Digit Med*, vol. 3, pp. 127, 2020.
- [35] R. Mc Ardle, S. Del Din, P. Donaghy, B. Galna, A. J. Thomas, and L. Rochester, "The Impact of Environment on Gait Assessment: Considerations from Real-World Gait Analysis in Dementia Subtypes," *Sensors*, 21, 2021].
- [36] S. Del Din, A. Godfrey, C. Mazzà, S. Lord, and L. Rochester, "Free-living monitoring of Parkinson's disease: Lessons from the field," *Mov Disord*, vol. 31, no. 9, pp. 1293-313, Sep, 2016.
- [37] M. D. Rinderknecht, M. Zanon, T. A. Boonstra, L. Angelini, D. Stanev, G. G. Chan, L. Bunn, F. Dondelinger, R. Hosking, J. Freeman, J. Hobart, J. Marsden, and L. Craveiro, "An observational study to assess validity and reliability of smartphone sensor-based gait and balance assessments in multiple sclerosis: Floodlight GaitLab protocol," *DIGITAL HEALTH*, vol. 9, pp. 20552076231205284, 2023/01/01, 2023.
- [38] C. H. Polman, S. C. Reingold, B. Banwell, M. Clanet, J. A. Cohen, M. Filippi, K. Fujihara, E. Havrdova, M. Hutchinson, L. Kappos, F. D. Lublin, X. Montalban, P. O'Connor, M. Sandberg-Wollheim, A. J. Thompson, E. Waubant, B. Weinshenker, and J. S. Wolinsky, "Diagnostic criteria for multiple sclerosis: 2010 Revisions to the McDonald criteria," *Annals of Neurology*, vol. 69, no. 2, pp. 292-302, 2011/02/01, 2011.
- [39] A. J. Thompson, B. L. Banwell, F. Barkhof, W. M. Carroll, T. Coetzee, G. Comi, J. Correale, F. Fazekas, M. Filippi, M. S. Freedman, K. Fujihara, S. L. Galetta, H. P. Hartung, L. Kappos, F. D. Lublin, R. A. Marrie, A. E. Miller, D. H. Miller, X. Montalban, E. M. Mowry, P. S. Sorensen, M. Tintore, A. L. Traboulsee, M. Trojano, B. M. J. Uitdehaag, S. Vukusic, E. Waubant, B. G. Weinshenker, S. C. Reingold, and J. A. Cohen, "Diagnosis

- of multiple sclerosis: 2017 revisions of the McDonald criteria," *Lancet Neurol*, vol. 17, no. 2, pp. 162-173, Feb, 2018.
- [40] E. Ibalá, S. Coupaud, and A. Kerr, "Comparison of the muscle pattern variability during treadmill walking (fixed and self-pace) and overground walking of able-bodied adults," *Journal of Annals of Bioengineering*, vol. 1, pp. 1-11, 2019.
  - [41] D. Stanev, R. Klimas, M. Płonka, N. Napiórkowski, L. Angelini, T. A. Boonstra, G. González Chan, J. Marsden, M. D. Rinderknecht, and M. Zanon, "Annotation and synchronization of smartphone inertial measurement unit and motion capture data collected in a laboratory setting to study gait and balance," *arXiv Preprint*, vol. 2411.06809, 2024.
  - [42] Motek, "HBM2 - Reference Manual. HBM-Gait. Available via: <https://knowledge.motekmedical.com/wp-content/uploads/2019/07/HBM2-Reference-Manual-Lower-body.pdf> [Accessed May 2025]," Motek, Amsterdam, Netherlands, 2018.
  - [43] L. Angelini, D. Y. Stanev, and M. Zanon. "Detection of Gait Activity. World Intellectual Property Organization (Patent No. WO/2024/170587). [https://patentscope.wipo.int/search/en/detail.jsf?docId=WO2024170587&\\_cid=P20-MAP55C-11648-1](https://patentscope.wipo.int/search/en/detail.jsf?docId=WO2024170587&_cid=P20-MAP55C-11648-1)."
  - [44] A. Godfrey, A. K. Bourke, G. M. O'laighin, P. van de Ven, and J. Nelson, "Activity classification using a single chest mounted tri-axial accelerometer," *Med Eng Phys*, vol. 33, no. 9, pp. 1127-35, Nov, 2011.
  - [45] S. O. Madgwick, A. J. Harrison, and A. Vaidyanathan, "Estimation of IMU and MARG orientation using a gradient descent algorithm," *IEEE Int Conf Rehabil Robot*, vol. 2011, pp. 5975346, 2011.
  - [46] F. Kluge, S. Del Din, A. Cereatti, H. Gaßner, C. Hansen, J. L. Helbostad, J. Klucken, A. Küderle, A. Müller, L. Rochester, M. Ullrich, B. M. Eskofier, and C. Mazzà, "Consensus based framework for digital mobility monitoring," *PLoS One*, vol. 16, no. 8, pp. e0256541, 2021.
  - [47] M. El-Gohary, S. Pearson, J. McNames, M. Mancini, F. Horak, S. Mellone, and L. Chiari, "Continuous monitoring of turning in patients with movement disability," *Sensors (Basel)*, vol. 14, no. 1, pp. 356-69, Dec 27, 2013.
  - [48] S. Bonnet, C. Bassompierre, C. Godin, S. Lesecq, and A. Barraud, "Calibration methods for inertial and magnetic sensors," *Sensors and Actuators A: Physical*, vol. 156, no. 2, pp. 302-311, 2009/12/01/, 2009.
  - [49] M. Nazarahari, and H. Rouhani, "Semi-Automatic Sensor-to-Body Calibration of Inertial Sensors on Lower Limb Using Gait Recording," *IEEE Sensors Journal*, vol. 19, no. 24, pp. 12465-12474, 2019.
  - [50] S. Wold, K. Esbensen, and P. Geladi, "Principal component analysis," *Chemometrics and Intelligent Laboratory Systems*, vol. 2, no. 1, pp. 37-52, 1987/08/01/, 1987.
  - [51] J. McCamley, M. Donati, E. Grimpampi, and C. Mazzà, "An enhanced estimate of initial contact and final contact instants of time using lower trunk inertial sensor data," *Gait Posture*, vol. 36, no. 2, pp. 316-8, Jun, 2012.
  - [52] G. R. Lee, R. Gommers, F. Waselewski, K. Wohlfahrt, and A. O'Leary, "PyWavelets: A Python package for wavelet analysis," *Journal of Open Source Software*, vol. 4, no. 36, pp. 1237, 2019.
  - [53] D. A. Winter, *Biomechanics and Motor Control of Human Movement*, 4th ed.: John Wiley & Sons, Inc., 2009.
  - [54] A. A. Taha, and A. Hanbury, "Metrics for evaluating 3D medical image segmentation: analysis, selection, and tool," *BMC Med Imaging*, vol. 15, pp. 29, Aug 12, 2015.
  - [55] R. McElreath, *Statistical Rethinking: A Bayesian Course with Examples in R and STAN*, 2nd ed.: Chapman and Hall/CRC, 2020.
  - [56] L. Rochester, C. Mazzà, A. Mueller, B. Caulfield, M. McCarthy, C. Becker, R. Miller, P. Piraino, M. Viceconti, W. P. Dartee, J. Garcia-Aymerich, A. A. Aydemir, B. Vereijken, V.

- Arnera, N. Ammour, M. Jackson, T. Hache, and R. Roubenoff, "A Roadmap to Inform Development, Validation and Approval of Digital Mobility Outcomes: The Mobilise-D Approach," *Digit Biomark*, vol. 4, no. Suppl 1, pp. 13-27, Winter, 2020.
- [57] A. Paraschiv-Ionescu, A. Soltani, and K. Aminian, "Real-world speed estimation using single trunk IMU: methodological challenges for impaired gait patterns()," *Annu Int Conf IEEE Eng Med Biol Soc*, vol. 2020, pp. 4596-4599, Jul, 2020.
- [58] A. Paraschiv-Ionescu, C. J. Newman, L. Carcreff, C. N. Gerber, S. Armand, and K. Aminian, "Locomotion and cadence detection using a single trunk-fixed accelerometer: validity for children with cerebral palsy in daily life-like conditions," *J Neuroeng Rehabil*, vol. 16, no. 1, pp. 24, Feb 4, 2019.
- [59] P. Abry, *Ondelettes et turbulences: multirésolutions, algorithmes de décomposition, invariance d'échelle et signaux de pression*, Paris, France: Diderot multimédia éd, 1997.
- [60] S. H. Shin, and C. G. Park, "Adaptive step length estimation algorithm using optimal parameters and movement status awareness," *Med Eng Phys*, vol. 33, no. 9, pp. 1064-71, Nov, 2011.
- [61] M. E. Mico-Amigo, T. Bonci, A. Paraschiv-Ionescu, M. Ullrich, C. Kirk, A. Soltani, A. Kuderle, E. Gazit, F. Salis, L. Alcock, K. Aminian, C. Becker, S. Bertuletti, P. Brown, E. Buckley, A. Cantu, A. E. Carsin, M. Caruso, B. Caulfield, A. Cereatti, L. Chiari, I. D'Ascanio, B. Eskofier, S. Fernstad, M. Froehlich, J. Garcia-Aymerich, C. Hansen, J. M. Hausdorff, H. Hiden, E. Hume, A. Keogh, F. Kluge, S. Koch, W. Maetzler, D. Megaritis, A. Mueller, M. Niessen, L. Palmerini, L. Schwickert, K. Scott, B. Sharrack, H. Sillen, D. Singleton, B. Vereijken, I. Vogiatzis, A. J. Yarnall, L. Rochester, C. Mazza, S. Del Din, and D. c. Mobilise, "Assessing real-world gait with digital technology? Validation, insights and recommendations from the Mobilise-D consortium," *J Neuroeng Rehabil*, vol. 20, no. 1, pp. 78, Jun 14, 2023.
- [62] A. Hickey, S. Del Din, L. Rochester, and A. Godfrey, "Detecting free-living steps and walking bouts: validating an algorithm for macro gait analysis," *Physiol Meas*, vol. 38, no. 1, pp. N1-n15, Jan, 2017.
- [63] N. Ji, H. Zhou, K. Guo, O. W. Samuel, Z. Huang, L. Xu, and G. Li, "Appropriate Mother Wavelets for Continuous Gait Event Detection Based on Time-Frequency Analysis for Hemiplegic and Healthy Individuals," *Sensors (Basel)*, vol. 19, no. 16, Aug 8, 2019.
- [64] M. Roerdink, B. H. Coolen, B. H. Clairbois, C. J. Lamoth, and P. J. Beek, "Online gait event detection using a large force platform embedded in a treadmill," *J Biomech*, vol. 41, no. 12, pp. 2628-32, Aug 28, 2008.
- [65] J. A. Zeni, Jr., J. G. Richards, and J. S. Higginson, "Two simple methods for determining gait events during treadmill and overground walking using kinematic data," *Gait Posture*, vol. 27, no. 4, pp. 710-4, May, 2008.
- [66] I. Härdi, S. A. Bridenbaugh, Y. J. Gschwind, and R. W. Kressig, "The effect of three different types of walking aids on spatio-temporal gait parameters in community-dwelling older adults," *Aging Clin Exp Res*, vol. 26, no. 2, pp. 221-8, Apr, 2014.
- [67] M. Mundt, J. P. Batista, B. Markert, C. Bollheimer, and T. Laurentius, "Walking with rollator: a systematic review of gait parameters in older persons," *Eur Rev Aging Phys Act*, vol. 16, pp. 15, 2019.

# Declarations

## Ethics approval and consent to participate

All participants provided written informed consent to participate in the study, which was approved by the UK Health Research Authority (IRAS Project ID: 302099).

## Data availability

Requesting the data underlying this publication requires a detailed, hypothesis-driven statistical analysis plan collaboratively developed by the requester and company subject matter experts.

Such requests should be directed to [dbm.datarequest@roche.com](mailto:dbm.datarequest@roche.com) for consideration.

Anonymized records for individual patients across >1 data source external to Roche cannot, and should not, be linked due to a potential increase in the risk of patient reidentification.

## Competing interests

The author(s) declared the following potential conflicts of interest concerning the research, authorship, and/or publication of this article: LA is an employee of F. Hoffmann-La Roche Ltd. DS is a consultant for F. Hoffmann-La Roche Ltd via Capgemini Engineering. MP is a contractor for Roche Polska Sp. z o.o. RK is a contractor for F. Hoffmann-La Roche Ltd. NN is a contractor for Roche Polska Sp. z o.o. GGC is an employee of the University of Plymouth, which received research funding from F. Hoffman-La Roche Ltd. LB is an employee of the University of Plymouth, which received research funding for this study from F. Hoffmann-La Roche Ltd. PSG is an employee of the University of Plymouth. RH has received indirect funding from F. Hoffmann-La Roche Ltd via a University of Plymouth Enterprise Ltd contract. JF is an employee of the University of Plymouth, which received research funding for this study from F. Hoffmann-

La Roche Ltd. JM is an employee of the University of Plymouth, which received research funding for this study from F. Hoffmann-La Roche Ltd. JH or affiliated institutions have received either consulting fees, honoraria, support to attend meetings, clinical service support, or research support from Acorda, Bayer Schering Pharma, Biogen Idec., Brickell Biotech, F. Hoffmann-La Roche Ltd, Global Blood Therapeutics, Sanofi Genzyme, Merck Serono, Novartis, Oxford Health Policy Forum, Teva, and Vantia. LC is an employee of, and shareholder in F. Hoffmann-La Roche Ltd. MDR is a contractor for F. Hoffmann-La Roche Ltd. MZ is an employee of and shareholder in F. Hoffmann-La Roche Ltd.

## Funding

This research was funded by F. Hoffmann-La Roche Ltd, Basel, Switzerland.

## Authors' contributions

LA: Conceptualization, Data Curation, Formal analysis, Methods, Software, Visualization, Writing - Review & Editing

DS: Conceptualization, Data Curation, Formal analysis, Methods, Software, Visualization, Writing - Review & Editing

MP: Data Curation, Formal analysis, Methods, Software, Visualization, Writing - Review & Editing

RK: Data Curation, Formal analysis, Methods, Software, Visualization, Writing - Review & Editing

NN: Data Curation, Formal analysis, Software, Visualization, Writing - Review & Editing

GGC: Data Curation, Investigation, Writing - Review & Editing

LB: Data Curation, Investigation, Writing - Review & Editing

PSG: Data Curation, Investigation, Writing - Review & Editing

RH: Data Curation, Investigation, Writing - Review & Editing

JF: Data Curation, Investigation, Writing - Review & Editing, Conceptualization

JM: Conceptualization, Data Curation, Investigation, Supervision, Writing - Review & Editing

JH: Data Curation, Investigation, Writing - Review & Editing, Conceptualization

LC: Conceptualization, Funding acquisition, Methods Writing - Review & Editing

MDR: Conceptualization, Formal analysis, Methods, Project administration, Supervision, Writing  
– Original Draft, Writing - Review & Editing

MZ: Conceptualization, Formal analysis, Methods, Project administration, Supervision, Writing -  
Review & Editing

## Acknowledgements

The authors would like to thank all study participants and their families. The authors would also like to thank the following current and former colleagues at F. Hoffmann-La Roche Ltd for their contributions to and support of the study: Pawel Biernacki, Matthias Bobst, Marco Bosatta, Sandro Fritz, Francisco Gavidia, Petra Hauser, Wiktoria Kasprzyk, Dominik Kedziora, Piotr Konorski, Kostas Kritsas, Hugo Le Gall, Dominik Lenart, Michael Lindemann, Florian Lipsmeier, David Meyer, Madalina Ogica, Emanuele Passerini, Szymon Piwowar, Elena Spyridou, and Lukasz Syzmankiewicz. The authors would also like to thank colleagues at the University of Plymouth for their contributions and support to the study: Michael Paisey, Tanya King, Joanne

Lind, Carol Lunn, and Maureen Pedder. Editorial and writing support was provided by Sven Holm, PhD, contractor for F. Hoffmann-La Roche Ltd.

## Supplementary appendix

### Data synchronisation

The work by Stanev et al. [1] presents a methodology for synchronising data across various measurement systems, enhancing the accuracy and reliability of data acquisition in clinical settings. Smartphones and Gait Up IMUs were synchronised using a mechanical perturbation method, where devices were subjected to a rocking motion to generate identifiable gyroscope signal patterns, allowing for temporal offset correction. A vibration-based method was also employed to synchronise smartphones and Gait Up IMUs with mocap system, using a vibrating motor to generate pulse sequences for alignment. Despite potential lags introduced by the vibrating motor's inertia or pulse sequence delays, post hoc synchronisation techniques, such as the force-based lag estimation method, were utilised to correct these discrepancies. This method leverages the relationship between ground reaction forces (GRFs) and center of mass acceleration during walking or standing, using accelerometer readings from a smartphone to estimate and rectify lags, thereby improving data quality.

### Gait event detection from the reference tools

The mocap and GRF data were processed using the Gait Offline Analysis Tool (GOAT) software package, which was developed by manufacturer [2]. This software offers built-in algorithms that can accurately identify initial and final contacts. GOAT features two distinct event detection algorithms, namely the kinematic-based and GRF-based algorithms, making it a reliable tool for gait analysis in healthy individuals and those with pathological conditions [3]. The kinematic-based algorithm, also known as the coordinate-based algorithm, detects gait events based on the positions of sacrum and foot markers filtered with a 2nd order Butterworth filter and cut-off frequency of 6Hz. This algorithm was validated against the events identified from the vertical GRF, which is considered the gold standard for such analysis. Comprehensive details of the algorithmic implementation have been previously described [4]. The GRF-based algorithm, on the other hand, relies on force data that have been converted to center-of-pressure (COP) data. This conversion

enables the identification of gait events using the characteristic COP “butterfly” pattern that emerges over time. The initial and final contacts obtained from the COP profile were validated against those based on kinematic data, which were assumed to represent “true” gait events [5]. By default, the gait events used in the subsequent analyses were identified from the gold standard using the GRF-based algorithm. However, when a walking aid was used, the kinematic-based method was employed to identify gait events.

The gait events were obtained from the silver standard using a developed gait analysis package, which has been previously validated against motion capture systems in both young and elderly individuals [6-9] and in limping walking conditions across a range of walking speeds between 0.9 and 2.0 m/second [10]. Foot-worn IMUs were primarily used to detect gait events (e.g., initial and final contacts). In contrast, the IMU attached to the lower back was used for more accurate turn detection.

## Bayesian model

The model's outcome variable was the F1 Score, which ranged from 0 to 100 or from 0 to 1 when normalised. Given the bounded nature of the outcome variable, a beta distribution between 0 and 1 was employed. The beta distribution, which is utilised for modelling the F1 score, is parameterised by the  $\alpha$  and  $\beta$  parameters, which are positive. The  $\alpha$  and  $\beta$  parameters are parametrised in terms of  $\mu$  and  $k$  as follows (also known as beta regression reparametrisation):


$$\alpha = \mu \cdot k$$

$$\beta = (1 - \mu) \cdot k \text{ with } 0 < \mu < 1$$

A mathematical description of the Bayesian beta regression model is presented in Supplementary Figure S2.

The parameter  $\mu$  was utilised to link the predictor variables to the distribution parameters. Given the bounded constraint of  $\mu$  ( $0 < \mu < 1$ ), a link sigmoid function (logit) was introduced, leading to the formulation of a generalised linear model. Parameters  $a$  and  $b$  represent the intercept and age parameters;  $s$ , and  $d$  denote the sex and disease parameters corresponding to their index variables;  $u$  accommodates subject intercepts to account for repeated measures within subjects;  $e$  and  $h$  represent the effects under examination (e.g., indoor/outdoor and walking aid). Furthermore, parameter  $d$  was constrained to behave as an ordinal predictor variable as a function of disease severity. Tables S1 and S2 summarise the contrast parameter estimates of the Bayesian model for the initial and final contacts, respectively.

## Supplementary figures and tables

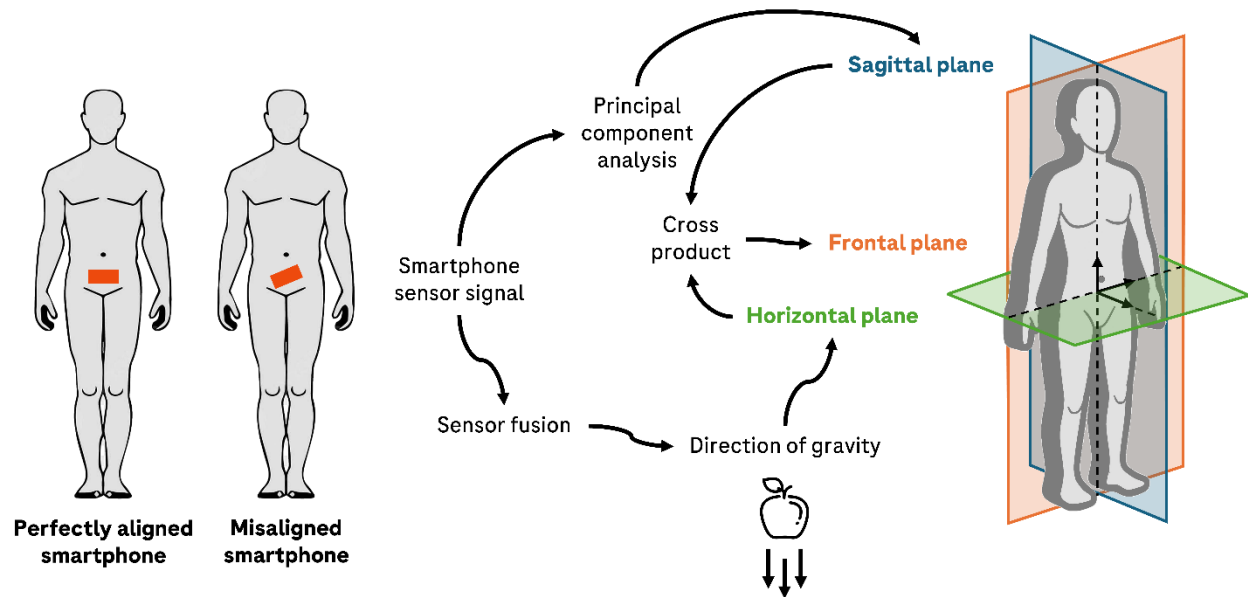

**Supplementary Figure S1. Alignment with gravity and anatomical reference frame.** Smartphone accelerometers that are not aligned with the direction of movement may negatively impact the performance of the step detection algorithm. Thus, aligning the smartphone accelerometer and gyroscope signal to the direction of gravity and the anatomical reference frame that is agnostic to the smartphone's orientation is required to ensure the algorithm's ability to handle differently oriented smartphone devices. The direction of gravity, or horizontal plane, was identified through sensor fusion (i.e., Madgwick filter [11] applied to the smartphone accelerometer and gyroscope signals). The two remaining planes that together with the horizontal plane make up the orientation-agnostic anatomical reference frame were identified as follows: The sagittal plane (i.e., the direction of movement) was identified by conducting principal component analysis on the gravity-aligned accelerometer data and the frontal plane through computing the cross product of the horizontal and sagittal planes.

$$\begin{aligned}
\mathbf{Fl}_{obs}[i] &\sim \text{Beta}(\alpha[i], \beta[i]), \quad i \in \{0, \dots, N-1\} \\
\alpha[i] &= \mu[i] \cdot \kappa \\
\beta[i] &= (1 - \mu[i]) \cdot \kappa \\
\text{logit}(\mu[i]) &= a + b \cdot \bar{age} + s[\text{sex}_{idx}[i]] + d \cdot \sum_{j=0}^{\text{disease}_{idx}[i]-1} \delta_j + u[\text{subject}_{idx}[i]] + e[\text{environment}_{idx}[i]] + h[\text{aid}_{idx}[i]] \\
\kappa &\sim \text{HalfCauchy}(20) \\
a &\sim \text{Normal}(1, 1) \\
s &= 1/100 \\
b &\sim \text{Normal}(0, s) \\
\mu_{sub} &\sim \text{Normal}(0, s) \\
\sigma_{sub} &\sim \text{HalfCauchy}(s) \\
s[j] &\sim \text{Normal}(0, s), \quad j \in \{\text{M}, \text{F}\} \\
d &\sim \text{Normal}(0, s) \\
\delta[j] &\sim \text{Diriclet}(1), \quad j \in \{\text{disease categories}\} \\
u[j] &\sim \text{Normal}(\mu_{sub}, \sigma_{sub}), \quad j \in \{\text{subjects}\} \\
e[j] &\sim \text{Normal}(0, s), \quad j \in \{\text{Indoor}, \text{Outdoor}\} \\
h[j] &\sim \text{Normal}(0, s), \quad j \in \{\text{W/Aid}, \text{W/O Aid}\}
\end{aligned}$$

**Supplementary Figure S2. Complete description of the Bayesian model.**

**Supplementary Table S1. Algorithm performance for initial contact.**

Cohort <sup>a</sup>	Smartphone location	Precision								Recall								F1 score							
		Mdn	P05	Q1	Q3	P95	ws-	Mean [95% CI]		Mdn	P05	Q1	Q3	P95	ws-	Mean [95% CI]		Mdn	P05	Q1	Q3	P95	ws-	Mean [95% CI]	
	IQR								IQR								IQR								
Structured in-lab treadmill 2MWT – smartphone vs mocap																									
HC	BB	100	99.0	100	100	100	0.0	99.8 [99.4 - 100]		100	100	100	100	100	0.0	100 [100 - 100]		100	99.5	100	100	100	0.0	99.9 [99.7 - 100]	
	BF	100	98.2	100	100	100	0.0	99.7 [99.4 - 100]		100	99.9	100	100	100	0.0	100 [100 - 100]		100	98.9	100	100	100	0.0	99.9 [99.7 - 100]	
	PB - left	100	98.49	100	100	100	0.0	99.7 [99.4 - 100]		100	99.5	100	100	100	0.0	99.8 [99.6 - 100]		100	98.0	100	100	100	0.0	99.8 [99.6 - 100]	
	PB - right	100	7.6	100	100	100	0.0	99.7 [99.3 - 100]		100	99.5	100	100	100	0.0	99.9 [99.8 - 100]		100	98.2	100	100	100	0.0	99.8 [99.6 - 100]	
	PF - left	100	98.7	100	100	100	0.0	99.7 [99.4 - 100]		100	100	100	100	100	0.0	100 [100 - 100]		100	99.2	100	100	100	0.0	99.9 [99.9 - 100]	
	PF - right	100	99.0	100	100	100	0.0	99.8 [99.4 - 100]		100	99.8	100	100	100	0.0	100 [100 - 100]		100	99.2	100	100	100	0.0	99.9 [99.9 - 100]	
Mild	BB	100	99.29	100	100	100	0.0	99.9 [99.6 - 100]		100	99.1	100	100	100	0.0	99.9 [99.7 - 100]		100	99.3	100	100	100	0.0	99.9 [99.7 - 100]	
PwMS	BF	100	9.6	100	100	100	0.0	99.8 [99.6 - 100]		100	99.5	100	100	100	0.0	99.9 [99.7 - 100]		100	99.7	100	100	100	0.0	99.9 [99.6 - 100]	
	PB - left	100	99.2	100	100	100	0.0	99.8 [99.5 - 100]		100	98.1	100	100	100	0.0	99.8 [99.5 - 100]		100	99.0	100	100	100	0.0	99.8 [99.5 - 100]	
	PB - right	100	99.6	100	100	100	0.0	99.9 [99.6 - 100]		100	98.8	100	100	100	0.0	99.8 [99.6 - 100]		100	99.3	100	100	100	0.0	99.8 [99.6 - 100]	
	PF - left	100	100	100	100	100	0.0	100 [99.9 - 100]		100	99.6	100	100	100	0.0	99.9 [99.8 - 100]		100	99.6	100	100	100	0.0	99.9 [99.9 - 100]	
	PF - right	100	99.6	100	100	100	0.0	100 [99.9 - 100]		100	100	100	100	100	0.0	100 [100 - 100]		100	99.8	100	100	100	0.0	100 [100 - 100]	
Moderate	BB	100	97.4	100	100	100	0.0	99.5 [98.9 - 100]		100	99.3	100	100	100	0.0	99.8 [99.6 - 100]		100	97.5	100	100	100	0.0	99.7 [99.3 - 100]	
PwMS	BF	100	94.4	100	100	100	0.0	99.3 [98.4 - 100]		100	98.1	100	100	100	0.0	99.6 [99.2 - 100]		100	95.6	99.9	100	100	0.1	99.4 [98.9 - 100]	
	PB - left	100	91.5	100	100	100	0.0	98.9 [97.9 - 100]		100	90.1	99.8	100	100	0.2	98.0 [95.9 - 100]		100	91.6	99.6	100	100	0.4	98.3 [96.9 - 99.8]	
	PB - right	100	92.7	100	100	100	0.0	99.1 [98.3 - 100]		100	93.9	100	100	100	0.0	99.2 [98.6 - 99.9]		100	94.8	99.8	100	100	0.2	99.2 [98.5 - 99.9]	
	PF - left	100	92.6	100	100	100	0.0	98.3 [95.8 - 100]		100	97.7	100	100	100	0.0	98.6 [96.2 - 100]		100	96.0	100	100	100	0.0	98.4 [96.0 - 100]	
	PF - right	100	97.0	100	100	100	0.0	99.4 [98.8 - 100]		100	98.8	100	100	100	0.0	99.8 [99.7 - 100]		100	97.9	99.9	100	100	0.1	99.6 [99.3 - 100]	

Severe	BB	100	88.7	97.4	100	100	2.6	94.4 [86.8 - 100]	100	85.8	97.4	100	100	2.6	94.3 [86.7 - 100]	99.7	88.8	97.4	100	100	2.6	94.4 [86.8 - 100]
PwMS	BF	100	90.8	99.3	100	100	0.7	97.3 [93.5 - 100]	100	97.8	99.4	100	100	0.6	99.6 [99.3 - 99.9]	100	94.1	99.5	100	100	0.5	98.2 [95.7 - 100]
	PB - left	99.8	87.5	97.4	100	100	2.6	95.6 [91.2 - 99.9]	100	86.8	95.7	100	100	4.3	95.5 [91.2 - 99.8]	99.7	87.6	96.9	100	100	3.1	95.5 [91.2 - 99.9]
	PB - right	100	78.8	98.4	100	100	1.6	96.3 [92.6 - 100]	100	78.1	97.9	100	100	2.1	95.5 [91.7 - 99.4]	99.9	79.0	98.0	100	100	2.0	95.9 [92.1 - 99.6]
	PF - left	100	65.2	98.9	100	100	1.1	95.3 [89.9 - 100]	100	66.4	98.5	100	100	1.5	95.5 [90.1 - 100]	100	65.8	98.8	100	100	1.2	95.4 [90.0 - 100]
	PF - right	100	70.3	98.7	100	100	1.3	95.9 [92.0 - 99.8]	100	69.8	97.7	100	100	2.3	95.8 [91.6 - 100]	99.9	70.2	98.8	100	100	1.2	95.8 [91.8 - 99.9]
All PwMS	BB	100	92.7	100	100	100	0.0	98.1 [95.9 - 100]	100	95.8	100	100	100	0.9	98.2 [96.0 - 100]	100	94.8	99.7	100	100	0.3	98.2 [95.9 - 100]
	BF	100	95.7	100	100	100	0.0	98.9 [97.7 - 100]	100	97.8	100	100	100	0.9	99.7 [99.5 - 99.9]	100	95.8	99.8	100	100	0.2	99.2 [98.4 - 100]
	PB - left	100	90.6	99.8	100	100	0.2	98.2 [96.8 - 99.6]	100	87.6	99.5	100	100	0.5	97.8 [96.2 - 99.3]	100	88.6	99.4	100	100	0.6	97.9 [96.5 - 99.4]
	PB - right	100	91.5	100	100	100	0.0	98.5 [97.3 - 99.7]	100	91.6	99.7	100	100	0.3	98.3 [97.1 - 99.5]	100	91.1	99.7	100	100	0.3	98.4 [97.2 - 99.6]
	PF - left	100	91.1	100	100	100	0.0	97.9 [96.0 - 99.8]	100	97.3	100	100	100	0.9	98.1 [96.2 - 100]	100	95.3	99.8	100	100	0.2	98.0 [96.1 - 99.8]
	PF - right	100	91.5	100	100	100	0.0	98.5 [97.3 - 99.9]	100	96.2	100	100	100	0.9	98.7 [97.4 - 100]	100	95.5	99.8	100	100	0.2	98.6 [97.4 - 99.9]
<i>Structured in-lab treadmill 2MWT – smartphone vs Gait Up</i>																						
HC	BB	100	100	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]
	BF	100	100	100	100	100	0.0	100 [99.9 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	99.9	100	100	100	0.0	100 [100 - 100]
	PB - left	100	100	100	100	100	0.0	100 [99.9 - 100]	100	99.5	100	100	100	0.0	99.8 [99.6 - 100]	100	99.5	100	100	100	0.0	99.9 [99.7 - 100]
	PB - right	100	100	100	100	100	0.0	99.9 [99.8 - 100]	100	99.7	100	100	100	0.0	100 [99.9 - 100]	100	99.9	100	100	100	0.0	100 [99.9 - 100]
	PF - left	100	99.6	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	99.7	100	100	100	0.0	100 [99.9 - 100]
	PF - right	100	100	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]
Mild	BB	100	99.1	100	100	100	0.0	99.8 [99.5 - 100]	100	99.0	100	100	100	0.0	99.8 [99.6 - 100]	100	99.5	100	100	100	0.0	99.8 [99.5 - 100]
PwMS	BF	100	99.6	100	100	100	0.0	99.8 [99.5 - 100]	100	99.4	100	100	100	0.0	99.9 [99.6 - 100]	100	99.7	100	100	100	0.0	99.8 [99.6 - 100]
	PB - left	100	99.6	100	100	100	0.0	99.8 [99.4 - 100]	100	98.0	100	100	100	0.0	99.8 [99.5 - 100]	100	99.0	100	100	100	0.0	99.8 [99.4 - 100]
	PB - right	100	99.6	100	100	100	0.0	99.8 [99.4 - 100]	100	98.7	100	100	100	0.0	99.7 [99.4 - 100]	100	99.3	100	100	100	0.0	99.7 [99.4 - 100]
	PF - left	100	99.6	100	100	100	0.0	99.9 [99.9 - 100]	100	99.7	100	100	100	0.0	100 [99.9 - 100]	100	99.6	100	100	100	0.0	99.9 [99.9 - 100]
	PF - right	100	99.6	100	100	100	0.0	99.9 [99.8 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	99.8	100	100	100	0.0	100 [99.9 - 100]

Moderate	BB	100	97.9	99.6	100	100	0.4	99.6	[99.3 - 99.9]	100	99.0	100	100	100	0.0	99.8	[99.7 - 100]	100	98.5	99.8	100	100	0.2	99.7	[99.5 - 99.9]
PwMS	BF	100	98.2	100	100	100	0.0	99.7	[99.5 - 99.9]	100	99.5	100	100	100	0.0	99.9	[99.7 - 100]	100	98.9	100	100	100	0.0	99.8	[99.6 - 99.9]
	PB - left	100	96.5	99.3	100	100	0.7	99.2	[98.7 - 99.8]	100	92.6	100	100	100	0.0	98.3	[96.1 - 100]	100	94.3	99.4	100	100	0.6	98.6	[97.3 - 100]
	PB - right	100	97.0	99.3	100	100	0.7	99.4	[99.0 - 99.8]	100	96.2	99.5	100	100	0.5	99.4	[98.9 - 99.8]	100	97.3	99.1	100	100	0.9	99.4	[99.0 - 99.8]
	PF - left	100	97.9	99.7	100	100	0.3	98.3	[95.5 - 100]	100	98.5	100	100	100	0.0	98.6	[96.1 - 100]	100	98.7	99.5	100	100	0.5	98.4	[95.8 - 100]
	PF - right	100	97.8	100	100	100	0.0	99.5	[99.2 - 99.9]	100	98.4	100	100	100	0.0	99.7	[99.6 - 99.9]	100	99.4	99.4	100	100	0.6	99.6	[99.4 - 99.9]
Severe	BB	100	93.2	97.8	100	100	2.2	94.6	[86.5 - 100]	100	86.2	100	100	100	1.6	94.2	[86.1 - 100]	99.8	89.6	97.7	100	100	2.3	94.4	[86.3 - 100]
PwMS	BF	100	92.4	99.5	100	100	0.5	97.5	[93.8 - 100]	100	99.5	100	100	100	0.5	99.5	[98.8 - 100]	100	92.6	99.5	100	100	0.5	98.3	[95.7 - 100]
	PB - left	99.8	84.1	98.4	100	100	1.6	95.6	[90.5 - 100]	99.8	97.7	100	100	100	2.3	95.2	[90.1 - 100]	99.5	83.5	98.2	100	100	1.8	95.4	[90.3 - 100]
	PB - right	100	90.2	99.0	100	100	1.0	97.6	[95.0 - 100]	100	99.0	100	100	100	1.0	96.6	[93.7 - 99.6]	100	87.8	99.0	100	100	1.0	97.1	[94.4 - 99.8]
	PF - left	100	67.8	99.0	100	100	1.0	96.3	[91.9 - 100]	100	98.9	100	100	100	1.1	96.3	[91.9 - 100]	100	68.6	99.2	100	100	0.8	96.3	[91.9 - 100]
	PF - right	100	84.8	99.3	100	100	0.7	97.8	[95.8 - 99.9]	100	98.7	100	100	100	1.3	97.4	[95.1 - 99.8]	100	84.6	98.7	100	100	1.3	97.6	[95.4 - 99.9]
All PwMS	BB	100	95.8	99.5	100	100	0.5	98.2	[95.7 - 100]	100	96.5	100	100	100	0.0	98.1	[95.6 - 100]	100	96.8	99.5	100	100	0.5	98.2	[95.7 - 100]
	BF	100	97.2	100	100	100	0.0	99.0	[97.9 - 100]	100	98.9	100	100	100	0.0	99.7	[99.5 - 100]	100	98.3	99.8	100	100	0.2	99.3	[98.5 - 100]
	PB - left	100	92.4	99.1	100	100	0.9	98.3	[96.7 - 99.9]	100	90.2	100	100	100	0.0	97.7	[95.9 - 99.5]	100	90.8	99.3	100	100	0.7	98.0	[96.3 - 99.6]
	PB - right	100	95.6	99.6	100	100	0.4	99.0	[98.1 - 99.8]	100	91.4	99.5	100	100	0.5	98.6	[97.7 - 99.6]	100	95.2	99.3	100	100	0.7	98.8	[97.9 - 99.7]
	PF - left	100	97.9	99.8	100	100	0.2	98.1	[96.3 - 99.9]	100	97.9	100	100	100	0.0	98.2	[96.5 - 100]	100	98.0	99.6	100	100	0.4	98.2	[96.4 - 100]
	PF - right	100	97.4	100	100	100	0.0	99.1	[98.4 - 99.8]	100	97.4	100	100	100	0.0	99.1	[98.4 - 99.9]	100	98.1	99.7	100	100	0.3	99.1	[98.4 - 99.8]
<i>Structured in-lab overground 2MWT – smartphone vs Gait Up</i>																									
HC	BB	99.5	95.9	98.5	100	100	1.5	98.9	[98.4 - 99.4]	99.4	97.2	98.5	99.9	100	1.4	99.0	[98.6 - 99.3]	99.1	97.4	98.4	99.7	100	1/3	98.9	[98.7 - 99.2]
	BF	99.4	96.2	97.4	100	100	2.6	98.9	[98.3 - 99.3]	99.5	97.5	98.9	100	100	1.1	99.3	[99.0 - 100]	99.2	97.5	98.4	99.7	100	1.3	99.0	[98.7 - 99.3]
	PB - left	98.9	92.2	97.0	100	100	3.0	98.0	[97.0 - 98.9]	99.2	95.5	98.3	99.5	100	1.2	98.6	[98.0 - 99.2]	98.7	95.2	97.4	99.4	98.8	2.0	98.3	[97.7 - 98.8]
	PB - right	98.9	93.0	96.9	99.9	100	3.0	97.9	[97.1 - 98.7]	99.4	97.0	98.5	99.9	100	1.4	98.9	[98.5 - 99.3]	98.7	96.0	98.2	99.4	99.9	1.2	98.4	[97.9 - 98.9]
	PF - left	98.3	91.4	95.8	99.5	100	3.7	97.3	[96.3 - 98.3]	98.7	96.8	97.9	99.5	100	1.6	98.7	[98.3 - 99.1]	98.2	94.5	97.2	99.3	99.8	2.2	98.0	[97.4 - 98.6]
	PF - right	99.0	93.3	97.2	100	100	2.8	98.0	[97.2 - 98.8]	99.8	96.5	98.9	100	100	1.1	99.2	[98.8 - 99.6]	99.1	96.3	97.9	99.7	100	1.8	98.6	[98.1 - 99.0]

Mild	BB	98.8	96.7	97.8	100	100	2.2	98.6	[98.0 - 99.3]	99.4	97.2	97.7	100	100	2.3	98.9	[98.4 - 99.4]	98.9	96.9	98.5	99.4	100	0.9	98.9	[98.3 - 99.2]
PwMS	BF	98.4	95.9	97.7	100	100	2.3	98.4	[97.7 - 99.1]	100	97.7	99.1	100	100	0.9	99.4	[99.0 - 99.8]	99.1	97.7	98.4	99.4	100	1.1	98.9	[98.4 - 99.4]
	PB - left	97.7	93.3	96.4	98.9	100	2.5	97.3	[96.3 - 98.4]	98.8	95.0	97.7	100	100	2.3	98.4	[97.7 - 99.1]	97.6	96.3	97.3	98.6	99.5	1.4	97.8	[97.3 - 98.4]
	PB - right	97.6	93.2	96.0	99.5	100	3.4	97.4	[96.4 - 98.5]	97.8	96.3	97.0	100	100	3.0	98.2	[97.5 - 98.9]	97.7	95.2	97.0	99.1	99.7	2.1	97.8	[97.1 - 98.5]
	PF - left	97.5	93.2	96.4	98.5	100	2.0	97.2	[96.2 - 98.3]	99.4	93.2	96.7	100	100	3.3	98.2	[97.2 - 99.2]	97.9	94.9	96.5	99.2	99.8	2.7	97.7	[96.9 - 98.5]
	PF - right	98.1	92.2	96.0	99.0	100	2.9	97.2	[96.1 - 98.3]	99.0	96.4	98.3	99.5	100	1.2	98.7	[98.3 - 99.2]	98.0	95.5	97.4	99.2	99.7	1.8	97.9	[97.3 - 98.5]
Moderate	BB	100	95.4	98.1	100	100	1.9	98.6	[97.6 - 99.7]	99.0	97.2	98.4	100	100	1.6	99.0	[98.6 - 99.3]	99.3	96.7	98.9	99.7	100	0.8	98.8	[98.1 - 99.5]
PwMS	BF	100	92.4	97.9	100	100	2.1	98.1	[96.6 - 99.7]	100	95.1	98.8	100	100	1.2	99.0	[98.3 - 99.7]	99.4	94.1	98.7	100	100	1.3	98.5	[97.5 - 99.6]
	PB - left	100	93.1	97.9	100	100	2.1	98.3	[97.2 - 99.4]	98.7	95.6	97.8	99.4	100	1.6	98.4	[97.9 - 98.8]	98.9	95.4	97.9	99.4	99.7	1.5	98.3	[97.6 - 99.0]
	PB - right	100	93.3	97.7	100	100	2.3	97.9	[96.3 - 99.6]	98.6	94.8	97.4	98.9	100	1.5	97.7	[96.5 - 98.9]	98.7	92.4	98.0	99.4	100	1.4	97.8	[96.5 - 99.0]
	PF - left	99.3	93.2	96.4	100	100	3.6	97.8	[96.5 - 99.1]	99.4	96.9	98.4	100	100	1.6	99.0	[98.7 - 99.4]	99.3	95.6	97.7	99.7	100	2.0	98.4	[97.6 - 99.2]
	PF - right	98.9	93.1	97.8	100	100	2.2	97.7	[96.2 - 99.1]	99.4	96.9	98.8	100	100	1.2	99.0	[98.5 - 99.5]	99.0	95.5	98.2	99.5	100	1.3	98.4	[97.4 - 99.3]
Severe	BB	99.0	96.0	98.2	100	100	1.8	96.7	[92.2 - 100]	99.3	97.8	98.6	100	100	1.4	97.0	[92.4 - 100]	99.2	97.6	98.5	99.6	100	1.1	96.8	[92.3 - 100]
PwMS	BF	98.7	96.9	98.4	100	100	1.6	97.6	[95.0 - 100]	99.4	97.8	99.2	100	100	0.8	98.3	[96.1 - 100]	99.2	97.9	98.9	99.7	100	0.8	98.0	[95.5 - 100]
	PB - left	98.7	96.6	97.7	100	100	2.3	98.4	[97.6 - 99.2]	98.6	95.9	98.1	100	100	1.9	98.6	[98.0 - 99.1]	98.8	97.2	98.0	99.3	100	1.4	98.5	[97.9 - 99.1]
	PB - right	98.4	93.1	97.7	100	100	2.3	98.0	[96.9 - 99.1]	98.5	92.4	97.5	99.1	100	1.7	97.7	[96.9 - 98.5]	98.6	93.5	98.0	99.1	99.7	1.1	97.8	[96.9 - 98.8]
	PF - left	98.6	95.0	96.9	99.4	100	2.5	98.1	[97.4 - 98.8]	99.4	94.5	98.6	100	100	1.4	98.7	[98.1 - 99.4]	98.8	95.1	97.8	99.7	100	1.9	98.4	[97.8 - 99.0]
	PF - right	98.6	87.5	97.0	100	100	3.0	97.2	[95.3 - 99.0]	98.8	90.6	97.6	99.4	100	1.8	97.6	[96.2 - 99.0]	98.5	88.8	97.4	99.6	100	2.2	97.3	[95.7 - 99.0]
All PwMS	BB	99.3	95.7	97.9	100	100	2.1	97.9	[96.2 - 99.6]	99.3	97.2	98.4	100	100	1.6	98.2	[96.5 - 99.9]	99.2	97.0	98.5	99.6	100	1.1	98.0	[96.4 - 99.7]
	BF	99.3	94.8	98.2	100	100	1.8	98.0	[96.9 - 99.2]	100	96.8	99.1	100	100	0.9	98.8	[98.0 - 99.7]	99.2	97.1	98.8	99.7	100	1.0	98.4	[97.4 - 99.4]
	PB - left	98.8	93.4	97.2	100	100	2.8	98.1	[97.5 - 98.6]	98.7	95.1	97.9	99.5	100	1.5	98.4	[98.1 - 98.8]	98.6	96.0	97.6	99.4	100	1.7	98.2	[97.9 - 98.6]
	PB - right	98.9	92.9	97.1	100	100	2.9	97.8	[97.1 - 98.6]	98.4	93.6	97.3	99.3	100	2.1	97.8	[97.2 - 98.4]	98.5	93.8	97.3	99.3	99.9	1.9	97.8	[97.2 - 98.4]
	PF - left	98.6	93.2	96.5	100	100	3.5	97.8	[97.2 - 98.4]	99.4	94.9	98.3	100	100	1.7	98.7	[98.3 - 99.1]	98.7	95.0	97.6	99.7	100	2.1	98.2	[97.8 - 98.6]
	PF - right	98.6	92.4	96.4	100	100	3.6	97.4	[96.4 - 98.3]	99.1	95.6	98.1	100	100	1.9	98.4	[97.8 - 99.0]	98.6	94.6	97.6	99.5	100	1.8	97.9	[97.1 - 98.6]

### Structured real-world walking – smartphone vs Gait Up

HC	BF suggested location	100	98.7	99.7	100	100	0.5	99.7	[99.6 - 99.9]	100	98.8	100	100	100	0.0	99.8	[99.7 - 100]	100	98.7	99.7	100	100	0.6	99.8	[99.6 - 99.9]
Mild PwMS	BF suggested location	100	97.4	100	100	100	0.4	99.6	[99.1 - 100]	100	98.1	100	100	100	0.0	99.2	[97.7 - 100]	100	97.4	99.9	100	100	0.3	99.4	[98.4 - 100]
Moderate PwMS	BF suggested location	99.6	93.8	97.9	100	100	1.9	98.3	[97.3 - 99.4]	99.9	93.3	98.4	100	100	1.2	98.7	[97.9 - 99.5]	99.4	94.1	98.2	100	100	1.8	98.5	[97.6 - 99.3]
Severe PwMS	BF suggested location	98.0	73.0	90.6	100	100	8.2	93.8	[90.3 - 97.3]	98.9	80.3	96.4	100	100	5.4	95.4	[92.5 - 98.4]	98.5	76.4	93.6	99.6	100	7.2	94.4	[91.2 - 97.6]
All PwMS	BF suggested location	99.8	84.9	97.5	100	100	3.1	96.9	[95.5, 98.4]	100	85.9	98.4	100	100	1.5	97.6	[96.4 - 98.9]	99.5	87.0	98.1	100	100	2.4	97.2	[95.8 - 98.5]

### Unstructured real-world walking– smartphone vs Gait Up

HC	User-preferred location	97.4	92.7	94.8	98.4	99.2	5.1	96.5	[95.6 - 97.4]	98.7	94.8	97.4	99.5	99.8	3.2	98.1	[97.4 - 98.8]	97.8	92.6	96.2	98.8	99.4	4.5	97.1	[96.3 - 97.9]
Mild PwMS	User-preferred location	95.2	87.0	91.2	97.7	99.6	5.8	94.0	[91.7 - 96.3]	97.8	94.3	95.5	98.5	99.8	3.8	97.3	[96.5 - 98.1]	96.0	88.8	94.0	97.6	99.7	4.4	95.3	[93.7 - 96.9]
Moderate PwMS	User-preferred location	93.4	80.6	91.7	95.4	99.0	6.5	92.1	[89.8 - 94.3]	96.9	76.0	94.2	98.4	99.8	5.8	94.4	[91.7 - 97.1]	94.8	76.1	92.6	96.8	99.1	5.9	93.0	[90.6 - 95.4]
Severe PwMS	User-preferred location	89.6	73.0	83.3	93.1	97.1	11.4	88.1	[85.5 - 90.8]	93.9	80.9	90.3	96.8	99.1	9.3	92.3	[90.0 - 94.6]	91.8	78.5	86.6	94.7	97.4	10.1	89.9	[87.6 - 92.2]
All PwMS	User-preferred location	93.0	77.5	88.9	95.9	98.8	8.0	91.1	[89.6 - 92.6]	96.4	80.9	93.4	98.4	99.7	6.5	94.3	[92.9 - 95.7]	94.4	79.0	90.4	96.6	99.0	7.1	92.4	[91.0 - 93.8]

<sup>a</sup>Mild PwMS are PwMS with mild EDDS scores (EDSS: ≤ 3.5), moderate PwMS are PwMS with moderate EDSS scores (EDSS: 4.0 - 5.5), and severe PwMS are PwMS with severe EDSS scores (EDSS: 5.5 - 6.0).

BB = belt back; BF = belt front; CI = confidence interval; EDSS = Expanded Disability Status Scale; IQR = within-subject interquartile range; Mdn = median; P05 = 5th percentile; P95 = 95th percentile; PB = pocket back; PF = pocket front; PwMS = people with multiple sclerosis; Q1 = first quantile; Q3 = third quantile; ws-IQR = within-subject interquartile range.

**Supplementary Table S2. Algorithm performance for final contacts.**

Cohort <sup>a</sup>	Smartphone location	Precision								Recall								F1 score							
		Mdn	P05	Q1	Q3	P95	ws-	Mean [95% CI]	Mdn	P05	Q1	Q3	P95	ws-	Mean [95% CI]	Mdn	P05	Q1	Q3	P95	ws-	Mean [95% CI]			
	IQR								IQR								IQR								
Structured in-lab treadmill 2MWT – smartphone vs mocap																									
HC	BB	100	99.0	100	100	100	0.0	99.8 [99.4 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	99.5	100	100	100	0.0	99.9 [99.7 - 100]			
	BF	100	99.0	100	100	100	0.0	99.8 [99.4 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	99.4	100	100	100	0.0	99.9 [99.7 - 100]			
	PB - left	100	98.7	100	100	100	0.0	99.8 [99.4 - 100]	100	99.2	100	100	100	0.0	99.8 [99.5 - 100]	100	98.0	100	100	100	0.0	99.8 [99.6 - 100]			
	PB - right	100	98.4	100	100	100	0.0	99.7 [99.4 - 100]	100	99.9	100	100	100	0.0	100 [99.9 - 100]	100	99.2	100	100	100	0.0	99.8 [99.7 - 100]			
	PF - left	100	99.0	100	100	100	0.0	99.8 [99.4 - 100]	100	99.9	100	100	100	0.0	100 [99.9 - 100]	100	99.2	100	100	100	0.0	99.9 [99.7 - 100]			
	PF - right	100	99.0	100	100	100	0.0	99.8 [99.4 - 100]	100	99.4	100	100	100	0.0	99.9 [99.8 - 100]	100	98.9	100	100	100	0.0	99.8 [99.7 - 100]			
Mild	BB	100	100	100	100	100	0.0	100 [99.9 - 100]	100	100	100	100	100	0.0	100 [99.9 - 100]	100	100	100	100	100	0.0	100 [99.9 - 100]			
PwMS	BF	100	99.6	100	100	100	0.0	100 [99.9 - 100]	100	99.3	100	100	100	0.0	99.9 [99.8 - 100]	100	99.7	100	100	100	0.0	99.9 [99.9 - 100]			
	PB - left	100	100	100	100	100	0.0	99.9 [99.8 - 100]	100	98.9	100	100	100	0.0	99.8 [99.7 - 100]	100	99.1	100	100	100	0.0	99.9 [99.8 - 100]			
	PB - right	100	100	100	100	100	0.0	99.9 [99.9 - 100]	100	99.3	100	100	100	0.0	99.9 [99.7 - 100]	100	99.1	100	100	100	0.0	99.9 [99.8 - 100]			
	PF - left	100	100	100	100	100	0.0	100 [99.9 - 100]	100	98.6	100	100	100	0.0	99.6 [99.1 - 100]	100	99.3	99.9	100	100	0.1	99.8 [99.5 - 100]			
	PF - right	100	100	100	100	100	0.0	100 [100 - 100]	100	99.7	100	100	100	0.0	100 [99.9 - 100]	100	99.8	100	100	100	0.0	100 [100 - 100]			
Moderate	BB	100	97.7	100	100	100	0.0	99.5 [99.0 - 100]	100	99.5	100	100	100	0.0	99.9 [99.6 - 100]	100	97.5	100	100	100	0.0	99.7 [99.4 - 100]			
PwMS	BF	100	97.3	100	100	100	0.0	99.5 [98.9 - 100]	100	98.7	100	100	100	0.0	99.7 [99.4 - 100]	100	97.5	100	100	100	0.0	99.6 [99.3 - 99.9]			
	PB - left	100	97.7	100	100	100	0.0	99.5 [98.9 - 100]	100	93.5	100	100	100	0.0	98.2 [96.1 - 100]	100	95.6	100	100	100	0.0	98.7 [97.4 - 99.8]			
	PB - right	100	97.1	100	100	100	0.0	99.5 [99.0 - 100]	100	94.7	99.6	100	100	0.4	99.3 [98.8 - 99.9]	100	96.7	99.8	100	100	0.2	99.4 [99.0 - 99.8]			
	PF - left	100	93.0	100	100	100	0.0	98.9 [97.5 - 100]	100	97.2	99.8	100	100	0.2	98.9 [97.5 - 100]	100	96.4	99.8	100	100	0.2	98.9 [97.5 - 100]			
	PF - right	100	97.3	100	100	100	0.0	99.5 [98.9 - 100]	100	98.9	100	100	100	0.0	99.8 [99.6 - 100]	100	98.0	100	100	100	0.0	99.6 [99.3 - 100]			

Severe	BB	100	90.6	99.3	100	100	0.6	94.7 [86.5 - 100]	100	80.2	99.0	100	100	1.0	94.2 [86.0 - 100]	100	85.2	99.0	100	100	1.0	94.4 [86.2 - 100]
PwMS	BF	100	61.8	99.0	100	100	1.0	93.2 [84.5 - 100]	100	95.3	99.4	100	100	0.6	95.3 [87.2 - 99.9]	99.9	75.7	99.0	100	100	1.0	94.0 [85.7 - 100]
	PB - left	100	13.4	98.8	100	100	1.2	90.4 [79.2 - 99.9]	100	12.6	98.5	100	100	1.5	90.1 [78.9 - 100]	99.7	13.0	98.3	100	100	1.7	90.2 [79.0 - 100]
	PB - right	100	95.9	99.1	100	100	0.9	99.1 [98.1 - 100]	100	90.5	98.9	100	100	1.1	97.8 [96.0 - 99.6]	99.8	92.0	98.9	100	100	1.1	98.4 [97.2 - 99.6]
	PF - left	100	92.6	99.1	100	100	0.9	97.7 [94.5 - 100]	99.6	92.6	99.3	100	100	0.7	96.8 [92.6 - 100]	99.7	95.2	98.7	100	100	1.3	97.2 [93.4 - 100]
	PF - right	100	92.5	99.2	100	100	0.8	98.7 [97.6 - 99.9]	100	90.4	97.5	100	100	2.5	98.0 [96.6 - 99.4]	99.5	93.5	98.2	100	100	1.8	98.4 [97.2 - 99.5]
All PwMS	BB	100	93.7	100	100	100	0.0	98.2 [95.8 - 100]	100	96.3	100	100	100	0.0	98.2 [95.8 - 100]	100	95.7	100	100	100	0.0	98.2 [95.8 - 100]
	BF	100	93.7	100	100	100	0.0	97.8 [95.2 - 100]	100	98.2	100	100	100	0.0	98.5 [96.1 - 100]	100	95.6	99.7	100	100	0.3	98.1 [95.6 - 100]
	PB - left	100	91.7	100	100	100	0.0	97.0 [93.6 - 100]	100	86.7	99.8	100	100	0.2	96.3 [92.8 - 99.8]	100	89.3	99.6	100	100	0.4	96.6 [93.2 - 100]
	PB - right	100	98.0	100	100	100	0.0	99.5 [99.1 - 99.9]	100	94.2	99.6	100	100	0.4	99.0 [98.4 - 99.7]	100	95.5	99.6	100	100	0.4	99.2 [98.8 - 99.7]
	PF - left	100	93.9	100	100	100	0.0	98.8 [97.7 - 100]	100	95.4	99.5	100	100	0.5	98.5 [97.1 - 99.9]	100	95.6	99.5	100	100	0.5	98.7 [97.4 - 99.9]
	PF - right	100	96.7	100	100	100	0.0	99.4 [99.0 - 99.8]	100	96.4	100	100	100	0.0	99.3 [98.9 - 99.8]	100	96.0	99.6	100	100	0.4	99.4 [99.0 - 99.7]

*Structured in-lab treadmill 2MWT – smartphone vs Gait Up*

HC	BB	100	100	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]
	BF	100	100	100	100	100	0.0	100 [99.9 - 100]	100	100	100	100	100	0.0	100 [100 - 100]	100	100	100	100	100	0.0	100 [100 - 100]
	PB - left	100	99.9	100	100	100	0.0	100 [99.9 - 100]	100	99.2	100	100	100	0.0	99.8 [99.5 - 100]	100	99.5	100	100	100	0.0	99.9 [99.7 - 100]
	PB - right	100	99.9	100	100	100	0.0	100 [99.9 - 100]	100	100	100	100	100	0.0	100 [99.9 - 100]	100	99.8	100	100	100	0.0	100 [99.9 - 100]
	PF - left	100	100	100	100	100	0.0	100 [99.9 - 100]	100	99.8	100	100	100	0.0	100 [99.9 - 100]	100	99.8	100	100	100	0.0	100 [100 - 100]
	PF - right	100	100	100	100	100	0.0	100 [100 - 100]	100	99.8	100	100	100	0.0	99.9 [99.9 - 100]	100	99.9	100	100	100	0.0	100 [99.9 - 100]
Mild	BB	100	99.6	100	100	100	0.0	99.9 [99.9 - 100]	100	99.6	100	100	100	0.0	99.9 [99.8 - 100]	100	99.6	100	100	100	0.0	99.9 [99.9 - 100]
PwMS	BF	100	99.6	100	100	100	0.0	99.9 [99.9 - 100]	100	99.3	100	100	100	0.0	99.9 [99.7 - 100]	100	99.6	100	100	100	0.0	99.9 [99.8 - 100]
	PB - left	100	99.6	100	100	100	0.0	99.9 [99.8 - 100]	100	98.8	100	100	100	0.0	99.8 [99.6 - 100]	100	99.0	100	100	100	0.0	99.9 [99.7 - 100]
	PB - right	100	99.6	100	100	100	0.0	99.9 [99.8 - 100]	100	99.2	100	100	100	0.0	99.8 [99.6 - 100]	100	99.0	100	100	100	0.0	99.9 [99.8 - 100]
	PF - left	100	99.6	100	100	100	0.0	99.9 [99.9 - 100]	100	99.0	100	100	100	0.0	99.6 [99.1 - 100]	100	99.5	100	100	100	0.0	99.8 [99.5 - 100]
	PF - right	100	99.6	100	100	100	0.0	100 [99.9 - 100]	100	99.7	100	100	100	0.0	99.9 [99.9 - 100]	100	99.8	100	100	100	0.0	100 [99.9 - 100]

Moderate	BB	100	97.3	100	100	100	0.0	99.5	[99.1 - 99.9]	100	98.7	100	100	100	0.0	99.8	[99.5 - 100]	100	97.3	99.8	100	100	0.2	99.6	[99.3 - 99.9]
PwMS	BF	100	96.3	99.6	100	100	0.4	99.3	[99.8 - 99.9]	100	98.0	100	100	100	0.0	99.7	[99.4 - 100]	100	96.4	99.8	100	100	0.2	99.5	[99.1 - 99.9]
	PB - left	100	96.4	99.6	100	100	0.4	99.3	[98.6 - 99.9]	100	92.9	99.6	100	100	0.4	98.2	[96.1 - 100]	100	95.9	99.2	100	100	0.8	98.7	[97.2 - 100]
	PB - right	100	97.3	100	100	100	0.0	99.5	[99.0 - 99.9]	100	94.6	99.6	100	100	0.4	99.2	[98.6 - 99.8]	100	96.7	99.5	100	100	0.5	99.3	[98.9 - 99.8]
	PF - left	100	97.4	100	100	100	0.0	98.3	[95.7 - 100]	100	97.0	99.7	100	100	0.3	98.4	[95.9 - 100]	100	97.9	99.6	100	100	0.4	98.4	[95.8 - 100]
	PF - right	100	97.6	100	100	100	0.0	99.5	[99.2 - 99.9]	100	98.3	100	100	100	0.0	99.7	[99.4 - 100]	100	98.4	99.6	100	100	0.4	99.6	[99.3 - 99.9]
Severe	BB	100	91.9	99	100	100	1.0	95.0	[86.9 - 100]	100	81.0	98.4	100	100	1.6	94.3	[86.1 - 100]	100	86.1	98.3	100	100	1.7	94.6	[86.5 - 100]
PwMS	BF	100	62.6	99.2	100	100	0.8	93.8	[85.4 - 100]	100	97.4	99.4	100	100	0.6	95.7	[87.9 - 100]	100	76.0	99.2	100	100	0.8	94.5	[86.5 - 100]
	PB - left	99.8	23.9	99.3	100	100	0.7	91.6	[81.3 - 100]	99.7	23.1	97.5	100	100	2.5	90.9	[80.6 - 100]	99.7	23.5	98.1	100	100	1.9	91.3	[81.0 - 100]
	PB - right	100	97.7	99.3	100	100	0.7	99.5	[99.1 - 99.8]	99.8	91.4	98.1	100	100	1.9	98.0	[96.3 - 99.6]	99.8	94.9	98.7	100	100	1.3	98.7	[97.7 - 99.6]
	PF - left	100	85.4	99.3	100	100	0.7	97.5	[95.0 - 99.9]	99.8	79.4	99.2	100	100	0.8	96.3	[92.4 - 100]	99.8	82.3	98.5	100	100	1.5	96.8	[93.6 - 100]
	PF - right	100	96.6	98.7	100	100	1.3	97.5	[93.7 - 99.9]	100	90.6	97.8	100	100	2.2	96.6	[92.6 - 100]	99.7	93.5	98.2	100	100	1.8	97.0	[93.1 - 99.9]
All PwMS	BB	100	97.1	99.8	100	100	0.2	98.2	[95.8 - 100]	100	96.2	100	100	100	0.0	98.1	[95.6 - 100]	100	97.2	99.7	100	100	0.3	98.2	[95.7 - 100]
	BF	100	96.9	99.6	100	100	0.4	97.8	[95.2 - 100]	100	98.6	100	100	100	0.0	98.5	[96.2 - 100]	100	96.6	99.7	100	100	0.3	98.1	[95.6 - 100]
	PB - left	100	95.2	99.8	100	100	0.2	97.1	[93.9 - 100]	100	90.8	99.3	100	100	0.7	96.4	[93.1 - 99.7]	100	94.2	99.2	100	100	0.8	96.7	[93.5 - 100]
	PB - right	100	97.6	99.8	100	100	0.2	99.6	[99.4 - 99.8]	100	94.1	99.5	100	100	0.5	99.0	[98.4 - 99.6]	100	96.4	99.5	100	100	0.5	99.3	[98.9 - 99.6]
	PF - left	100	95.6	99.8	100	100	0.2	98.5	[97.2 - 99.9]	100	93.2	99.5	100	100	0.5	98.1	[96.5 - 99.7]	100	95.8	99.5	100	100	0.5	98.3	[96.8 - 99.7]
	PF - right	100	97.1	100	100	100	0.0	99.0	[97.9 - 100]	100	96.4	99.6	100	100	0.4	98.8	[97.6 - 100]	100	97.5	99.6	100	100	0.4	98.9	[97.7 - 100]

*Structured in-lab overground 2MWT – smartphone vs Gait Up*

HC	BB	98.3	91.5	94.7	100	100	5.3	97.0	[96.0 - 98.1]	98.7	95.6	97.4	99.4	100	2.0	98.3	[97.8 - 98.8]	97.9	94.0	96.5	98.9	100	2.3	97.6	[97.0 - 98.2]
	BF	98.4	91.2	95.3	99.5	100	4.2	97.1	[96.1 - 98.2]	98.1	94.7	96.4	99.4	100	3.0	97.8	[97.2 - 98.4]	97.8	93.4	96.7	99.1	99.7	2.4	97.4	[96.8 - 98.1]
	PB - left	95.8	90.9	93.2	99.1	100	6.0	96.0	[94.9 - 97.1]	98.8	96.2	98.2	99.9	100	1.7	98.6	[98.1 - 99.2]	97.5	94.6	95.7	98.9	99.5	3.2	97.2	[96.6 - 97.9]
	PB - right	95.5	90.5	93.4	99.4	100	6.0	95.8	[94.6 - 97.0]	99.0	95.5	97.5	99.5	100	2.0	98.6	[98.1 - 99.0]	96.7	94.7	95.8	99.1	99.7	3.3	97.1	[96.5 - 97.8]
	PF - left	94.9	89.8	93.2	99.4	100	6.2	95.6	[94.5 - 96.7]	98.3	95.3	96.7	99.5	100	2.9	98.1	[97.5 - 98.7]	96.5	93.9	95.7	98.0	99.6	2.4	96.8	[96.2 - 97.4]
	PF - right	95.5	91.4	93.6	99.4	100	5.9	95.9	[94.8 - 97.0]	98.4	93.6	96.5	99.3	100	2.8	97.7	[96.9 - 98.5]	96.9	94.2	95.6	98.2	99.3	2.6	96.8	[96.1 - 97.5]

Mild	BB	96.2	93.0	94.9	99.4	100	4.5	96.6	[95.5 - 97.7]	98.2	96.4	97.5	98.9	100	1.3	98.2	[97.6 - 98.6]	97.2	95.0	96.4	98.5	99.7	2.1	97.3	[96.7 - 98.0]
PwMS	BF	96.1	92.9	94.7	99.4	100	4.6	96.7	[95.5 - 97.8]	98.0	95.0	96.0	99.4	100	3.4	97.8	[97.0 - 98.6]	97.2	94.3	96.3	98.7	99.7	2.4	97.2	[96.4 - 98.0]
	PB - left	95.0	91.6	93.2	96.2	100	3.0	95.2	[94.0 - 96.4]	99.4	95.6	97.8	100	100	2.2	98.7	[98.1 - 99.3]	97.0	95.3	96.2	97.9	98.8	1.7	96.9	[96.2 - 97.5]
	PB - right	95.2	91.6	93.1	97.2	100	4.1	95.3	[94.1 - 96.6]	98.9	96.3	98.4	100	100	1.6	98.8	[98.2 - 99.3]	96.7	94.8	95.9	98.0	99.7	2.1	97.0	[96.3 - 97.7]
	PF - left	94.5	91.4	93.1	96.7	100	3.6	94.9	[93.7 - 96.2]	98.8	91.4	95.9	99.4	100	3.5	97.6	[96.4 - 98.8]	96.4	91.7	95.1	97.7	99.3	2.6	96.2	[95.2 - 97.3]
	PF - right	95.3	91.1	93.8	97.6	100	3.8	95.4	[94.1 - 96.8]	98.4	93.4	97.5	99.1	100	1.6	98.0	[97.2 - 98.9]	96.8	94.6	96.1	97.7	99.5	1.6	96.7	[95.8 - 97.6]
Moderate	BB	99.2	92.4	96.0	100	100	4.0	97.5	[96.3 - 98.7]	98.8	94.5	97.6	99.4	100	1.8	98.2	[97.5 - 98.8]	98.1	95.5	97.3	99.1	99.7	1.8	97.8	[97.1 - 98.5]
PwMS	BF	99.2	90.9	95.6	100	100	4.4	97.1	[95.6 - 98.6]	98.5	93.2	97.0	99.4	100	2.4	97.9	[97.1 - 98.6]	98.0	93.2	97.0	99.0	99.6	2.0	97.4	[96.5 - 98.3]
	PB - left	98.5	93.0	95.5	100	100	4.5	97.4	[96.4 - 98.4]	98.5	94.6	97.3	99.4	100	2.1	98.0	[97.4 - 98.6]	98.0	95.5	96.8	98.8	99.4	2.0	97.7	[97.1 - 98.2]
	PB - right	98.7	92.4	95.4	99.4	100	4.0	97.0	[95.5 - 98.4]	98.3	94.1	97.1	99.3	100	2.2	97.2	[95.6 - 98.9]	97.8	90.7	97.1	98.9	99.5	1.8	97.0	[95.7 - 98.3]
	PF - left	98.2	92.4	95.4	98.9	100	3.6	96.8	[95.5 - 98.0]	98.6	94.3	97.2	99.4	100	2.2	98.1	[97.4 - 98.7]	98.2	94.6	96.7	98.8	99.5	2.1	97.4	[96.5 - 98.2]
	PF - right	98.4	92.3	95.2	99.4	100	4.3	96.9	[95.6 - 98.1]	99.0	94.7	97.2	99.5	100	2.2	98.4	[97.8 - 99.0]	97.7	95.5	96.6	99.2	99.7	2.6	97.6	[96.8 - 98.4]
Severe	BB	99.0	91.7	97.4	100	100	2.6	95.6	[90.3 - 100]	98.7	94.5	98.0	99.3	100	1.4	95.8	[90.5 - 100]	98.7	93.7	98.0	99.4	99.8	1.4	95.7	[90.4 - 100]
PwMS	BF	99.3	94.8	97.8	100	100	2.2	97.5	[95.2 - 99.9]	99.2	96.2	98.2	100	100	1.8	97.9	[95.9 - 99.8]	99.1	96.7	97.8	99.5	100	1.6	97.7	[95.6 - 99.8]
	PB - left	98.9	94.4	96.9	100	100	3.1	98.2	[97.5 - 99.0]	98.5	94.8	97.4	99.3	100	1.9	97.9	[97.2 - 98.6]	98.2	94.5	97.2	99.3	100	2.1	98.0	[97.4 - 98.6]
	PB - right	99.1	94.5	97.7	100	100	2.3	97.9	[96.6 - 99.2]	97.9	91.3	96.9	98.7	100	1.9	97.1	[95.9 - 98.3]	98.0	93.8	97.2	99.1	99.7	1.9	97.5	[96.4 - 98.6]
	PF - left	98.3	94.2	96.8	99.5	100	2.7	97.5	[96.5 - 98.6]	98.8	92.5	97.0	99.4	100	2.4	97.1	[95.1 - 99.1]	98.1	95.4	97.0	99.0	99.7	2.1	97.3	[95.8 - 98.7]
	PF - right	98.4	87.7	96.7	100	100	3.3	96.8	[94.9 - 98.7]	97.8	87.6	96.9	98.6	100	1.7	96.5	[94.5 - 98.5]	97.9	86.8	97.0	98.8	99.2	1.7	96.6	[94.7 - 98.5]
All PwMS	BB	98.7	92.8	96.1	100	100	3.9	96.6	[94.5 - 98.6]	98.6	95.0	97.7	99.3	100	1.7	97.3	[95.3 - 99.3]	98.2	95.1	97.1	99.2	99.7	2.1	96.9	[94.9 - 98.9]
	BF	98.6	92.8	95.9	100	100	4.1	97.2	[96.1 - 98.2]	98.7	94.4	97.2	99.4	100	2.2	97.8	[97.0 - 98.7]	98.4	94.3	97.0	99.2	99.9	2.2	97.5	[96.6 - 98.3]
	PB - left	98.1	92.5	95.2	100	100	4.8	97.2	[96.6 - 97.8]	98.6	94.9	97.6	99.4	100	1.8	98.1	[97.7 - 98.5]	97.9	94.6	96.7	98.8	99.7	2.2	97.6	[97.2 - 98.0]
	PB - right	98.4	92.0	95.4	99.4	100	4.0	96.9	[96.1 - 97.8]	98.4	95.1	97.1	99.4	100	2.2	97.5	[96.8 - 98.3]	97.7	93.9	96.7	98.9	99.7	2.2	97.2	[96.5 - 97.9]
	PF - left	97.4	92.0	94.7	99.3	100	4.6	96.6	[95.9 - 97.3]	98.8	92.4	96.9	99.4	100	2.6	97.6	[96.8 - 98.4]	97.7	93.9	96.2	98.7	99.7	2.4	97.1	[96.4 - 97.7]
	PF - right	97.6	91.3	95.2	99.4	100	4.2	96.5	[95.6 - 97.4]	98.4	94.0	97.2	99.4	100	2.2	97.6	[96.8 - 98.4]	97.7	94.7	96.6	98.9	99.7	2.3	97.0	[96.2 - 97.8]

Structured real-world walking – smartphone vs Gait Up																									
HC	BF suggested location	100	98.9	100	100	100	0.4	99.8	[99.6 - 100]	100	98.3	100	100	100	0.4	99.7	[99.5 - 100]	100	98.6	99.8	100	100	0.5	99.8	[99.6 - 99.9]
Mild PwMS	BF suggested location	100	96.2	100	100	100	0.4	99.4	[98.8 - 100]	100	96.8	100	100	100	0.0	99.1	[97.7 - 100]	100	96.5	99.9	100	100	0.3	99.0	[97.6 - 100]
Moderate PwMS	BF suggested location	99.5	93.4	98.2	100	100	1.7	98.3	[97.3 - 99.4]	99.5	92.3	97.3	100	100	1.8	98.3	[97.3 - 99.2]	99.4	94.2	97.5	100	100	2.1	98.3	[97.4 - 99.2]
Severe PwMS	BF suggested location	98.6	75.4	94.1	100	100	6.4	94.3	[90.7 - 97.9]	98.3	78.4	94.0	99.6	100	5.6	94.7	[91.5 - 97.9]	98.1	76.9	93.1	99.6	100	6.2	94.4	[91.0 - 97.8]
All PwMS	BF suggested location	99.8	86.4	97.4	100	100	2.8	97.1	[95.6 - 98.6]	99.5	87.4	97.4	100	100	2.3	97.1	[95.74 -98.5]	99.4	86.5	97.5	100	100	2.7	97.0	[95.6 - 98.4]
Unstructured real-world walking – smartphone vs Gait Up																									
HC	User-preferred location	97.5	92.8	95.5	98.5	99.3	4.8	96.7	[95.8 - 97.6]	98.8	92.9	97.1	99.3	99.6	3.8	97.7	[96.9 - 98.5]	97.8	92.6	95.8	98.6	99.5	4.7	97.0	[96.2 - 97.8]
Mild PwMS	User-preferred location	95.6	86.5	91.7	98.0	99.5	5.8	94.1	[91.7 - 96.6]	97.6	92.0	94.1	98.1	99.8	4.5	96.3	[95.1 - 97.5]	96.5	89.2	93.4	97.5	99.6	3.9	94.9	[93.1 - 96.8]
Moderate PwMS	User-preferred location	93.8	81.5	90.2	96.0	98.9	6.8	92.4	[90.3 - 94.5]	95.4	79.9	92.7	98.0	99.5	7.2	93.2	[90.6 - 95.8]	94.3	79.6	90.8	96.1	99.1	6.2	92.6	[90.4 - 94.9]
Severe PwMS	User-preferred location	90.0	75.2	86.5	94.2	97.9	9.9	89.1	[86.6 - 91.6]	92.6	78.1	87.5	97.2	98.7	11.8	91.3	[88.9 - 93.7]	90.9	79.4	85.5	94.6	97.2	9.7	89.8	[87.7 - 91.9]
All PwMS	User-preferred location	93.6	77.5	88.9	96.1	98.9	8.3	91.6	[90.2 - 93.0]	95.4	79.5	90.9	97.8	99.4	7.7	93.2	[91.8 - 94.7]	94.0	81.3	90.0	96.5	99.1	7.0	92.1	[90.8 - 93.5]

<sup>a</sup>Mild PwMS are PwMS with mild EDDS scores (EDSS: ≤ 3.5), moderate PwMS are PwMS with moderate EDSS scores (EDSS: 4.0 - 5.5), and severe PwMS are PwMS with severe EDSS scores (EDSS: 5.5 - 6.0).  
BB = belt back; BF = belt front; CI = confidence interval; EDSS = Expanded Disability Status Scale; IQR = within-subject interquartile range; Mdn = median; P05 = 5th percentile; P95 = 95th percentile; PB = pocket back; PF = pocket front; PwMS = people with multiple sclerosis; Q1 = first quantile; Q3 = third quantile; ws-IQR = within-subject interquartile range.

**Supplementary Table S3. Algorithm error measures for initial contacts.**

Cohort	Smartphone location	Constant error, mean [95% CI]	Absolute error, mean [95% CI]	Variable error, mean [95% CI]	Total variability, mean [95% CI]	Median error, mean [95% CI]	Median absolute error, mean [95% CI]	IQR error, mean [95% CI]
<i>Structured in-lab treadmill 2MWT – smartphone vs mocap</i>								
HC	BB	0.06 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.01 [0.01 - 0.02]	0.06 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.02 [0.02 - 0.02]
	BF	0.04 [0.03 - 0.05]	0.04 [0.03 - 0.05]	0.02 [0.01 - 0.02]	0.04 [0.04 - 0.05]	0.04 [0.03 - 0.05]	0.04 [0.03 - 0.05]	0.02 [0.02 - 0.03]
	PB - left	0.08 [0.07 - 0.09]	0.08 [0.07 - 0.09]	0.03 [0.03 - 0.03]	0.08 [0.08 - 0.09]	0.08 [0.07 - 0.09]	0.08 [0.07 - 0.09]	0.05 [0.04 - 0.06]
	PB - right	0.08 [0.07 - 0.08]	0.08 [0.07 - 0.08]	0.03 [0.03 - 0.04]	0.08 [0.08 - 0.09]	0.07 [0.07 - 0.08]	0.07 [0.07 - 0.08]	0.06 [0.06 - 0.07]
	PF - left	0.05 [0.04 - 0.05]	0.05 [0.04 - 0.06]	0.04 [0.03 - 0.04]	0.06 [0.05 - 0.07]	0.05 [0.04 - 0.05]	0.05 [0.04 - 0.05]	0.06 [0.06 - 0.07]
	PF - right	0.04 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.03 [0.03 - 0.03]	0.05 [0.05 - 0.06]	0.04 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.05 [0.05 - 0.06]
PwMS	BB	0.07 [0.06 - 0.08]	0.07 [0.07 - 0.08]	0.02 [0.02 - 0.02]	0.08 [0.07 - 0.08]	0.07 [0.06 - 0.08]	0.07 [0.06 - 0.08]	0.03 [0.02 - 0.03]
	BF	0.05 [0.04 - 0.06]	0.05 [0.05 - 0.06]	0.02 [0.02 - 0.03]	0.06 [0.05 - 0.07]	0.05 [0.04 - 0.06]	0.05 [0.04 - 0.06]	0.03 [0.03 - 0.03]
	PB - left	0.07 [0.06 - 0.09]	0.08 [0.08 - 0.09]	0.03 [0.03 - 0.03]	0.09 [0.08 - 0.1]	0.07 [0.06 - 0.08]	0.08 [0.07 - 0.09]	0.05 [0.04 - 0.05]
	PB - right	0.08 [0.07 - 0.08]	0.08 [0.07 - 0.09]	0.04 [0.03 - 0.04]	0.09 [0.08 - 0.09]	0.08 [0.07 - 0.08]	0.08 [0.07 - 0.08]	0.05 [0.05 - 0.06]
	PF - left	0.06 [0.05 - 0.07]	0.07 [0.06 - 0.08]	0.04 [0.03 - 0.04]	0.08 [0.07 - 0.08]	0.06 [0.05 - 0.07]	0.06 [0.06 - 0.07]	0.06 [0.05 - 0.07]
	PF - right	0.06 [0.05 - 0.06]	0.06 [0.05 - 0.07]	0.03 [0.03 - 0.04]	0.07 [0.06 - 0.07]	0.05 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.05 [0.05 - 0.06]

<i>Structured in-lab treadmill 2MWT – smartphone vs Gait Up</i>								
HC	BB	0.04 [0.03 - 0.05]	0.04 [0.03 - 0.05]	0.02 [0.01 - 0.02]	0.04 [0.04 - 0.05]	0.04 [0.03 - 0.05]	0.04 [0.03 - 0.05]	0.02 [0.02 - 0.02]
	BF	0.02 [0.01 - 0.03]	0.03 [0.02 - 0.03]	0.02 [0.02 - 0.02]	0.03 [0.02 - 0.04]	0.02 [0.01 - 0.03]	0.03 [0.02 - 0.03]	0.03 [0.02 - 0.03]
	PB - left	0.06 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.03 [0.03 - 0.03]	0.07 [0.06 - 0.08]	0.06 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.05 [0.04 - 0.06]
	PB - right	0.06 [0.05 - 0.07]	0.06 [0.06 - 0.07]	0.04 [0.03 - 0.04]	0.07 [0.07 - 0.08]	0.06 [0.05 - 0.07]	0.06 [0.05 - 0.07]	0.07 [0.06 - 0.07]
	PF - left	0.03 [0.02 - 0.04]	0.04 [0.04 - 0.05]	0.04 [0.03 - 0.04]	0.05 [0.04 - 0.05]	0.03 [0.02 - 0.03]	0.04 [0.03 - 0.04]	0.06 [0.06 - 0.07]
	PF - right	0.03 [0.02 - 0.03]	0.04 [0.03 - 0.04]	0.03 [0.03 - 0.03]	0.04 [0.04 - 0.05]	0.02 [0.02 - 0.03]	0.03 [0.03 - 0.04]	0.06 [0.05 - 0.06]
PwMS	BB	0.05 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.02 [0.02 - 0.03]	0.06 [0.05 - 0.07]	0.05 [0.05 - 0.06]	0.05 [0.05 - 0.06]	0.03 [0.03 - 0.03]
	BF	0.03 [0.02 - 0.04]	0.04 [0.04 - 0.05]	0.03 [0.02 - 0.03]	0.05 [0.04 - 0.06]	0.03 [0.02 - 0.04]	0.04 [0.03 - 0.05]	0.03 [0.03 - 0.04]
	PB - left	0.06 [0.04 - 0.07]	0.07 [0.06 - 0.08]	0.03 [0.03 - 0.04]	0.08 [0.07 - 0.09]	0.06 [0.04 - 0.07]	0.07 [0.06 - 0.08]	0.05 [0.04 - 0.05]
	PB - right	0.06 [0.06 - 0.07]	0.07 [0.06 - 0.07]	0.04 [0.03 - 0.04]	0.07 [0.07 - 0.08]	0.06 [0.05 - 0.07]	0.06 [0.06 - 0.07]	0.06 [0.05 - 0.06]
	PF - left	0.05 [0.04 - 0.06]	0.06 [0.05 - 0.07]	0.04 [0.03 - 0.04]	0.07 [0.06 - 0.07]	0.05 [0.04 - 0.06]	0.05 [0.05 - 0.06]	0.06 [0.05 - 0.07]
	PF - right	0.04 [0.03 - 0.05]	0.05 [0.04 - 0.06]	0.03 [0.03 - 0.04]	0.06 [0.05 - 0.06]	0.04 [0.03 - 0.05]	0.05 [0.04 - 0.05]	0.05 [0.05 - 0.06]
<i>Structured in-lab overground 2MWT – smartphone vs Gait Up</i>								
HC	BB	0.03 [0.03 - 0.04]	0.04 [0.03 - 0.04]	0.02 [0.02 - 0.02]	0.04 [0.03 - 0.04]	0.03 [0.03 - 0.04]	0.03 [0.03 - 0.04]	0.03 [0.02 - 0.03]
	BF	0.02 [0.01 - 0.02]	0.03 [0.02 - 0.03]	0.02 [0.02 - 0.03]	0.03 [0.03 - 0.04]	0.02 [0.01 - 0.02]	0.02 [0.02 - 0.03]	0.03 [0.02 - 0.04]
	PB - left	0.06 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.03 [0.03 - 0.04]	0.07 [0.06 - 0.07]	0.06 [0.05 - 0.07]	0.06 [0.05 - 0.07]	0.05 [0.05 - 0.06]
	PB - right	0.05 [0.04 - 0.06]	0.06 [0.05 - 0.06]	0.04 [0.04 - 0.04]	0.07 [0.06 - 0.07]	0.05 [0.04 - 0.06]	0.06 [0.05 - 0.06]	0.07 [0.06 - 0.08]

PwMS	PF - left	0.03 [0.03 - 0.04]	0.05 [0.05 - 0.05]	0.05 [0.04 - 0.05]	0.06 [0.05 - 0.06]	0.03 [0.02 - 0.04]	0.05 [0.04 - 0.05]	0.08 [0.08 - 0.09]
	PF - right	0.02 [0.02 - 0.03]	0.04 [0.04 - 0.04]	0.04 [0.04 - 0.04]	0.05 [0.04 - 0.05]	0.02 [0.01 - 0.03]	0.03 [0.03 - 0.04]	0.07 [0.06 - 0.07]
	BB	0.03 [0.03 - 0.04]	0.04 [0.04 - 0.04]	0.03 [0.02 - 0.03]	0.05 [0.04 - 0.05]	0.03 [0.03 - 0.04]	0.04 [0.03 - 0.04]	0.04 [0.03 - 0.04]
	BF	0.02 [0.02 - 0.03]	0.03 [0.03 - 0.04]	0.03 [0.02 - 0.03]	0.04 [0.04 - 0.04]	0.02 [0.02 - 0.03]	0.03 [0.03 - 0.04]	0.04 [0.03 - 0.04]
	PB - left	0.05 [0.04 - 0.05]	0.06 [0.05 - 0.06]	0.04 [0.04 - 0.04]	0.07 [0.06 - 0.07]	0.05 [0.04 - 0.05]	0.05 [0.05 - 0.06]	0.06 [0.05 - 0.06]
	PB - right	0.05 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.04 [0.04 - 0.05]	0.07 [0.07 - 0.07]	0.05 [0.04 - 0.05]	0.05 [0.05 - 0.06]	0.07 [0.06 - 0.07]
	PF - left	0.04 [0.03 - 0.04]	0.05 [0.05 - 0.05]	0.05 [0.04 - 0.05]	0.06 [0.06 - 0.07]	0.04 [0.03 - 0.04]	0.05 [0.04 - 0.05]	0.07 [0.07 - 0.08]
	PF - right	0.03 [0.02 - 0.03]	0.04 [0.04 - 0.05]	0.04 [0.04 - 0.04]	0.05 [0.05 - 0.06]	0.03 [0.02 - 0.03]	0.04 [0.03 - 0.04]	0.06 [0.06 - 0.07]
<i>Structured real-world walking – smartphone vs Gait Up</i>								
HC	BF suggested location	0.05 [0.05 - 0.06]	0.06 [0.05 - 0.06]	0.02 [0.02 - 0.03]	0.06 [0.06 - 0.07]	0.05 [0.05 - 0.06]	0.05 [0.05 - 0.06]	0.03 [0.03 - 0.04]
PwMS	BF suggested location	0.06 [0.05 - 0.06]	0.07 [0.06 - 0.07]	0.04 [0.03 - 0.04]	0.07 [0.07 - 0.08]	0.06 [0.05 - 0.06]	0.06 [0.06 - 0.07]	0.05 [0.04 - 0.05]
<i>Unstructured real-world walking – smartphone vs Gait Up</i>								
HC	User-preferred location	0.06 [0.05 - 0.06]	0.07 [0.06 - 0.08]	0.05 [0.04 - 0.05]	0.08 [0.07 - 0.09]	0.06 [0.05 - 0.06]	0.07 [0.06 - 0.07]	0.06 [0.05 - 0.08]
PwMS	User-preferred location	0.07 [0.06 - 0.07]	0.08 [0.08 - 0.09]	0.06 [0.06 - 0.06]	0.09 [0.09 - 0.1]	0.07 [0.06 - 0.07]	0.08 [0.07 - 0.08]	0.07 [0.07 - 0.08]

BB = belt back; BF = belt front; CI = confidence interval; IQR = interquartile range; PB = pocket back; PF = pocket front; PwMS = people with multiple sclerosis.

**Supplementary Table S4. Algorithm error measures for final contacts.**

Cohort	Smartphone location	Constant error, mean [95% CI]	Absolute error, mean [95% CI]	Variable error, mean [95% CI]	Total variability, mean [95% CI]	Median error, mean [95% CI]	Median absolute error, mean [95% CI]	IQR error, mean [95% CI]
<i>Structured in-lab treadmill 2MWT – smartphone vs mocap</i>								
HC	BB	0.03 [0.03 - 0.04]	0.03 [0.03 - 0.04]	0.02 [0.01 - 0.02]	0.04 [0.03 - 0.04]	0.03 [0.03 - 0.04]	0.03 [0.03 - 0.04]	0.02 [0.02 - 0.02]
	BF	0.02 [0.01 - 0.03]	0.03 [0.02 - 0.03]	0.02 [0.02 - 0.02]	0.03 [0.03 - 0.04]	0.02 [0.01 - 0.03]	0.03 [0.02 - 0.03]	0.03 [0.02 - 0.03]
	PB - left	0.04 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.02 [0.02 - 0.02]	0.05 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.03 [0.02 - 0.03]
	PB - right	0.04 [0.03 - 0.05]	0.04 [0.04 - 0.05]	0.03 [0.02 - 0.03]	0.05 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.05 [0.04 - 0.05]	0.04 [0.04 - 0.05]
	PF - left	0.02 [0.01 - 0.02]	0.04 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.05 [0.04 - 0.05]	0.03 [0.02 - 0.03]	0.04 [0.04 - 0.05]	0.08 [0.07 - 0.08]
	PF - right	0.01 [0.01 - 0.02]	0.04 [0.03 - 0.04]	0.04 [0.03 - 0.04]	0.04 [0.04 - 0.05]	0.02 [0.01 - 0.03]	0.04 [0.03 - 0.04]	0.07 [0.06 - 0.07]
PwMS	BB	0.01 [0.0 - 0.02]	0.04 [0.04 - 0.04]	0.03 [0.02 - 0.03]	0.04 [0.04 - 0.05]	0.01 [0.0 - 0.02]	0.04 [0.03 - 0.04]	0.03 [0.03 - 0.04]
	BF	-0.01 [-0.02 - -0.0]	0.04 [0.04 - 0.05]	0.03 [0.03 - 0.03]	0.05 [0.05 - 0.06]	-0.01 [-0.02 to -0.00]	0.04 [0.04 - 0.05]	0.04 [0.04 - 0.04]
	PB - left	0.02 [0.01 - 0.03]	0.05 [0.04 - 0.05]	0.03 [0.03 - 0.03]	0.05 [0.05 - 0.06]	0.02 [0.01 - 0.03]	0.05 [0.04 - 0.05]	0.04 [0.04 - 0.05]
	PB - right	0.01 [0.0 - 0.02]	0.05 [0.04 - 0.05]	0.04 [0.03 - 0.04]	0.05 [0.05 - 0.06]	0.01 [0.01 - 0.02]	0.05 [0.04 - 0.05]	0.06 [0.05 - 0.06]
	PF - left	0.01 [0.0 - 0.02]	0.05 [0.04 - 0.06]	0.05 [0.04 - 0.05]	0.06 [0.05 - 0.06]	0.02 [0.01 - 0.03]	0.05 [0.04 - 0.05]	0.07 [0.07 - 0.08]
	PF - right	0.00 [-0.01 - 0.0]	0.05 [0.04 - 0.05]	0.04 [0.04 - 0.05]	0.05 [0.05 - 0.06]	0.00 [-0.01 - 0.01]	0.04 [0.04 - 0.05]	0.07 [0.06 - 0.08]

*Structured in-lab treadmill 2MWT – smartphone vs Gait Up*

HC	BB	0.05 [0.04 - 0.06]	0.05 [0.04 - 0.06]	0.02 [0.02 - 0.02]	0.06 [0.05 - 0.06]	0.05 [0.04 - 0.06]	0.05 [0.04 - 0.06]	0.03 [0.02 - 0.03]
	BF	0.06 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.02 [0.02 - 0.02]	0.07 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.02 [0.02 - 0.03]
	PB - left	0.07 [0.07 - 0.08]	0.07 [0.07 - 0.08]	0.02 [0.02 - 0.02]	0.08 [0.07 - 0.08]	0.07 [0.07 - 0.08]	0.07 [0.07 - 0.08]	0.03 [0.03 - 0.04]
	PB - right	0.07 [0.07 - 0.08]	0.07 [0.07 - 0.08]	0.03 [0.02 - 0.03]	0.08 [0.07 - 0.08]	0.08 [0.07 - 0.08]	0.08 [0.07 - 0.08]	0.04 [0.04 - 0.05]
	PF - left	0.05 [0.04 - 0.05]	0.06 [0.05 - 0.06]	0.04 [0.04 - 0.05]	0.07 [0.06 - 0.07]	0.06 [0.05 - 0.06]	0.06 [0.06 - 0.07]	0.08 [0.07 - 0.09]
	PF - right	0.04 [0.04 - 0.05]	0.05 [0.05 - 0.06]	0.04 [0.04 - 0.04]	0.06 [0.06 - 0.07]	0.05 [0.04 - 0.06]	0.06 [0.05 - 0.06]	0.07 [0.06 - 0.08]
PwMS	BB	0.01 [0.0 - 0.03]	0.06 [0.05 - 0.06]	0.03 [0.03 - 0.04]	0.06 [0.06 - 0.07]	0.01 [0.0 - 0.03]	0.05 [0.05 - 0.06]	0.04 [0.04 - 0.05]
	BF	0.04 [0.03 - 0.05]	0.06 [0.05 - 0.07]	0.03 [nan - nan]	0.07 [0.06 - 0.07]	0.04 [0.03 - 0.05]	0.06 [0.05 - 0.07]	0.04 [0.04 - 0.04]
	PB - left	0.04 [0.02 - 0.05]	0.07 [0.06 - 0.08]	0.03 [0.03 - 0.04]	0.07 [0.07 - 0.08]	0.04 [0.03 - 0.06]	0.07 [0.06 - 0.08]	0.05 [0.04 - 0.05]
	PB - right	0.04 [0.03 - 0.05]	0.06 [0.06 - 0.07]	0.04 [0.03 - 0.04]	0.07 [0.06 - 0.08]	0.05 [0.03 - 0.06]	0.06 [0.05 - 0.07]	0.06 [0.05 - 0.06]
	PF - left	0.04 [0.03 - 0.05]	0.06 [0.05 - 0.07]	0.05 [0.04 - 0.05]	0.07 [0.06 - 0.08]	0.05 [0.04 - 0.06]	0.06 [0.05 - 0.07]	0.08 [0.07 - 0.08]
	PF - right	0.03 [0.02 - 0.03]	0.05 [0.05 - 0.06]	0.04 [0.04 - 0.05]	0.06 [0.06 - 0.07]	0.03 [0.02 - 0.04]	0.05 [0.05 - 0.06]	0.07 [0.06 - 0.07]

*Structured in-lab overground 2MWT – smartphone vs Gait Up*

HC	BB	0.06 [0.05 - 0.07]	0.06 [0.05 - 0.07]	0.02 [0.02 - 0.03]	0.06 [0.06 - 0.07]	0.06 [0.05 - 0.07]	0.06 [0.05 - 0.07]	0.03 [0.02 - 0.03]
	BF	0.06 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.02 [0.02 - 0.02]	0.07 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.06 [0.06 - 0.07]	0.03 [0.02 - 0.03]
	PB - left	0.07 [0.07 - 0.08]	0.07 [0.07 - 0.08]	0.03 [0.02 - 0.03]	0.08 [0.07 - 0.09]	0.07 [0.07 - 0.08]	0.07 [0.07 - 0.08]	0.04 [0.03 - 0.04]
	PB - right	0.07 [0.06 - 0.08]	0.07 [0.06 - 0.08]	0.03 [0.03 - 0.03]	0.08 [0.07 - 0.08]	0.07 [0.06 - 0.08]	0.07 [0.06 - 0.08]	0.05 [0.04 - 0.05]

PwMS	PF - left	0.06 [0.05 - 0.06]	0.06 [0.06 - 0.07]	0.04 [0.04 - 0.05]	0.07 [0.07 - 0.08]	0.06 [0.05 - 0.07]	0.06 [0.06 - 0.07]	0.08 [0.07 - 0.08]
	PF - right	0.05 [0.04 - 0.06]	0.05 [0.05 - 0.06]	0.04 [0.04 - 0.04]	0.06 [0.06 - 0.07]	0.05 [0.04 - 0.06]	0.05 [0.05 - 0.06]	0.07 [0.06 - 0.08]
	BB	0.02 [0.01 - 0.03]	0.05 [0.05 - 0.05]	0.03 [0.03 - 0.04]	0.06 [0.05 - 0.06]	0.02 [0.01 - 0.03]	0.05 [0.04 - 0.05]	0.05 [0.04 - 0.05]
	BF	0.03 [0.03 - 0.04]	0.05 [0.04 - 0.05]	0.03 [0.03 - 0.03]	0.06 [0.05 - 0.06]	0.03 [0.03 - 0.04]	0.05 [0.04 - 0.05]	0.04 [0.03 - 0.04]
	PB - left	0.04 [0.03 - 0.05]	0.06 [0.05 - 0.06]	0.04 [0.03 - 0.04]	0.06 [0.06 - 0.07]	0.04 [0.03 - 0.05]	0.05 [0.05 - 0.06]	0.05 [0.04 - 0.05]
	PB - right	0.04 [0.03 - 0.05]	0.06 [0.05 - 0.06]	0.04 [0.03 - 0.04]	0.07 [0.06 - 0.07]	0.04 [0.03 - 0.05]	0.05 [0.05 - 0.06]	0.05 [0.05 - 0.06]
	PF - left	0.03 [0.03 - 0.04]	0.06 [0.06 - 0.06]	0.05 [0.05 - 0.06]	0.07 [0.06 - 0.07]	0.03 [0.03 - 0.04]	0.06 [0.05 - 0.06]	0.09 [0.08 - 0.1]
	PF - right	0.03 [0.02 - 0.03]	0.05 [0.05 - 0.05]	0.04 [0.04 - 0.05]	0.06 [0.05 - 0.06]	0.03 [0.02 - 0.04]	0.05 [0.04 - 0.05]	0.07 [0.06 - 0.07]

#### *Structured real-world walking – smartphone vs Gait Up*

HC	BF suggested location	0.08 [0.08 - 0.09]	0.08 [0.08 - 0.09]	0.03 [0.02 - 0.03]	0.09 [0.08 - 0.1]	0.08 [0.08 - 0.09]	0.08 [0.08 - 0.09]	0.03 [0.03 - 0.04]
PwMS	BF suggested location	0.05 [0.04 - 0.06]	0.08 [0.07 - 0.08]	0.04 [0.04 - 0.05]	0.08 [0.08 - 0.09]	0.05 [0.04 - 0.07]	0.07 [0.07 - 0.08]	0.06 [0.05 - 0.06]

#### *Unstructured real-world walking – smartphone vs Gait Up*

HC	User-preferred location	0.07 [0.06 - 0.07]	0.08 [0.08 - 0.09]	0.06 [0.05 - 0.06]	0.09 [0.09 - 0.1]	0.07 [0.06 - 0.08]	0.08 [0.07 - 0.09]	0.07 [0.06 - 0.08]
PwMS	User-preferred location	0.05 [0.05 - 0.06]	0.08 [0.08 - 0.09]	0.07 [0.06 - 0.07]	0.1 [0.09 - 0.1]	0.06 [0.05 - 0.07]	0.08 [0.08 - 0.09]	0.09 [0.08 - 0.09]

BB = belt back; BF = belt front; CI = confidence interval; IQR = interquartile range; PB = pocket back; PF = pocket front; PwMS = people with multiple sclerosis.

**Supplementary Table S5. Parameter estimates of the Bayesian model examining the influence of different factors on the F1 score for initial contacts.**

<b>Parameter</b>	<b>Mean</b>	<b>Median</b>	<b>Std</b>	<b>5%</b>	<b>95%</b>	<b>IQR</b>	<b>Z- score</b>	<b>P &gt;  z </b>
Female - Male	-0.197	-0.204	1.393	-2.421	2.076	1.891	-0.142	0.887
Indoors - Outdoors	-1.788	-1.779	1.359	-4.021	0.463	1.839	-1.315	0.188
With walking aid vs without walking aid	-0.697	-0.719	1.409	-3.047	1.582	1.877	-0.495	0.621
Disease	-1.769	-1.77	1.01	-3.42	-0.123	1.373	-1.751	0.08

**Supplementary Table S6. Parameter estimates of the Bayesian model examining the influence of different factors on the F1 score for final contacts.**

<b>Parameter</b>	<b>Mean</b>	<b>Median</b>	<b>Std</b>	<b>5%</b>	<b>95%</b>	<b>IQR</b>	<b>Z-score</b>	<b>P &gt;  z </b>
Female vs male	-0.168	-0.174	1.388	-2.437	2.135	1.861	-0.121	0.904
Indoors vs outdoors	-2.061	-2.052	1.356	-4.308	0.169	1.834	-1.52	0.129
With walking aid vs without walking aid	-0.426	-0.421	1.419	-2.738	1.937	1.847	-0.3	0.764
Disease	-1.784	-1.791	1.0	-3.401	-0.171	1.369	-1.783	0.075

## References

- [1] D. Stanev, R. Klimas, M. Płonka, N. Napiórkowski, L. Angelini, T. A. Boonstra, G. González Chan, J. Marsden, M. D. Rinderknecht, and M. Zanon, "Annotation and synchronization of smartphone inertial measurement unit and motion capture data collected in a laboratory setting to study gait and balance," *arXiv Preprint*, vol. 2411.06809, 2024.
- [2] Motek, "HBM2 - Reference Manual. HBM-Gait. Available via: <https://knowledge.motekmedical.com/wp-content/uploads/2019/07/HBM2-Reference-Manual-Lower-body.pdf> [Accessed May 2025]," Motek, Amsterdam, Netherlands, 2018.
- [3] L. Aarts, S. Papegaaij, F. Steenbrink, and P. Martens, "Quality of Treadmill Embedded Force Plates for Gait Analysis. A White Paper. Available via: <https://knowledge.motekmedical.com/wp-content/uploads/2019/04/Whitepaper-ForcePlates-June2018.pdf> [Accessed May 2025]," Motek, Amsterdam, Netherlands, 2018.
- [4] J. A. Zeni, Jr., J. G. Richards, and J. S. Higginson, "Two simple methods for determining gait events during treadmill and overground walking using kinematic data," *Gait Posture*, vol. 27, no. 4, pp. 710-4, May, 2008.
- [5] M. Roerdink, B. H. Coolen, B. H. Clairbois, C. J. Lamoth, and P. J. Beek, "Online gait event detection using a large force platform embedded in a treadmill," *J Biomech*, vol. 41, no. 12, pp. 2628-32, Aug 28, 2008.
- [6] B. Mariani, C. Hoskovec, S. Rochat, C. Büla, J. Penders, and K. Aminian, "3D gait assessment in young and elderly subjects using foot-worn inertial sensors," *J Biomech*, vol. 43, no. 15, pp. 2999-3006, Nov 16, 2010.
- [7] B. Mariani, S. Rochat, C. J. Büla, and K. Aminian, "Heel and toe clearance estimation for gait analysis using wireless inertial sensors," *IEEE Trans Biomed Eng*, vol. 59, no. 11, pp. 3162-8, Nov, 2012.
- [8] B. Mariani, H. Rouhani, X. Crevoisier, and K. Aminian, "Quantitative estimation of foot-flat and stance phase of gait using foot-worn inertial sensors," *Gait Posture*, vol. 37, no. 2, pp. 229-34, Feb, 2013.
- [9] A. Brégou Bourgeois, B. Mariani, K. Aminian, P. Y. Zambelli, and C. J. Newman, "Spatio-temporal gait analysis in children with cerebral palsy using, foot-worn inertial sensors," *Gait Posture*, vol. 39, no. 1, pp. 436-42, 2014.
- [10] H. Schwameder, M. Andress, E. Graf, and G. Strutzenberger, "Validation of an Imu-System (Gait-Up) to Identify Gait Parameters in Normal and Induced Limping Walking Conditions." pp. 1130-1333.
- [11] S. O. Madgwick, A. J. Harrison, and A. Vaidyanathan, "Estimation of IMU and MARG orientation using a gradient descent algorithm," *IEEE Int Conf Rehabil Robot*, vol. 2011, pp. 5975346, 2011.