

qLOOK: A Minimal Information System for Digital Storage and Reproducible Analysis of qPCR experiments

Mirco Castoldi

Department of Gastroenterology, Hepatology and Infectious Diseases, Medical Faculty and University Hospital, Heinrich Heine University Düsseldorf, Germany

Center for Integrated Oncology Aachen Bonn Cologne Düsseldorf (CIO ABCD), Düsseldorf, Germany

Corresponding author:

Mirco Castoldi, mirco.castoldi@hhu.de

Keywords: *qLOOK, qPCR, Electronic laboratory notebooks, minimal information, MIQE*

Word count: *1846*

Number of figures and tables: *3*

Financial support: *This work was supported by a grant of the German Research Foundation (DFG, Deutsche Forschungsgemeinschaft) to M.C. (Reference: CA 830/3-1)*

Conflicts of interest: *The author has no relevant financial or non-financial interests to disclose*

Abstract

Objective: *Quantitative real-time PCR is widely used for gene expression analysis, yet inconsistencies in data storage and reporting limit reproducibility. While MIQE guidelines define the minimal information required for publication, they do not specify structured digital storage formats compatible with long-term reanalysis. This work presents qLOOK (qPCR-LOG-boOK), a tool for standardized digital storage and reproducible analysis of qPCR experiments.*

Results: *qLOOK is a modular R-based system that extracts data from Thermo Fisher/ABI .EDS files, formats it into a structured table (qLOOK_Data.xlsx), performs normalization and statistical analysis, and generates a log file (qLOOK_Summary.txt) recording reference genes, calibrators, and analytical parameters. All required R libraries are automatically installed and loaded, allowing researchers without coding experience to use the scripts. By preserving the qLOOK_Data table and the qLOOK_Summary log, users can reproduce or extend analyses without reprocessing raw files. While currently limited to .EDS files, the modular design allows adaptation to additional qPCR formats in the future.*

Besides providing an easy and transparent approach to analyze qPCR experiments, qLOOK also provides a minimal, standardized, and transparent solution for digital documentation, enhancing reproducibility, supporting long-term data stewardship, and facilitating integration into electronic laboratory notebooks or publication supplementary material.

Introduction

Quantitative real-time PCR (qPCR) is a cornerstone technique in molecular biology, enabling precise quantification of gene expression and nucleic acid abundance across diverse biological applications. Despite its widespread use, qPCR remains one of the most variable and inconsistently reported experimental methods. To improve transparency and reproducibility, the Minimum Information for Publication of Quantitative Real-Time PCR Experiments (MIQE) guidelines were introduced in 2009 (1, 2), providing standards for experimental design, data analysis, and reporting. While MIQE defines what information should be reported, it does not specify how these data should be digitally stored and preserved. In practice, qPCR data are often kept as instrument-specific output files, accompanied by manually curated spreadsheets or text documents. This unstructured approach poses several challenges: essential metadata are frequently missing or inconsistently recorded; raw data may only be accessible through proprietary software; and analytical steps such as normalization or calibrator are rarely documented in a reproducible way. Consequently, even within the same laboratory, reproducing or reanalyzing experiments conducted years earlier can be difficult.

The adoption of electronic laboratory notebooks (ELNs) and data management mandates has increased the need for standardized digital storage. ELNs provide a framework for traceable data capture but currently lack templates for qPCR data organization and processing. To address this gap, this work introduces qLOOK (qPCR-LOG-boOK), a minimal information system for digital storage, reproducible analysis, and reporting of qPCR experiments. qLOOK translates MIQE principles into a practical, script-based digital format compatible with ELNs and open-source environments. qLOOK consists of three modular R scripts that (i) extract essential experimental data from raw .EDS files (Thermo Fisher/ABI), (ii) perform normalization and statistical analysis, and (iii) generate visualizations and summaries. The core outputs are two structured text files: qLOOK_Data, containing formatted tabular qPCR data, and qLOOK_Summary, a log of all analytical steps. Together, these files form a reproducible,

platform-independent data package, allowing any user with access to R and the scripts to reproduce or extend the analysis, independent of proprietary software or hardware.

Beyond technical reproducibility, qLOOK supports open and digital research workflows. Including qLOOK outputs in ELNs or as supplementary material enhances transparency, facilitates data quality assessment, and aligns with FAIR (3) (Findable, Accessible, Interoperable, and Reusable) principles. The system is intentionally minimal, focusing on essential elements for reproducible storage and reanalysis. While currently supporting Thermo Fisher .EDS files, its structure can be readily adapted to additional formats, enabling broader standardization of qPCR data in digital research environments.

Design principles and scopae

qLOOK was designed to provide a minimal, reproducible, and broadly accessible system for analysis of qPCR experiments and digital storage. Its primary goal is to allow researchers to preserve essential experimental data and analytical steps in a structured format that enables full reanalysis, independent of proprietary software or instrument-specific tools. While compatibility with ELNs and FAIR data principles emerged as a beneficial feature, the core design focuses on simplicity, portability, and reproducibility.

The script works with EDS file generated by Thermo Fisher/ABI software and has been successfully tested on multiple cyclers, including 7500, 7500Fast, StepOnePlus, Viiia7, QuantStudio6 and QuantStudio3. R was selected as the implementation language to ensure platform independence and accessibility. The scripts are fully modular and require no coding experience; user input is handled via pop-up windows, and a future option for self-contained Windows executables is envisioned for users who do not wish to install R or RStudio, or working in environments with restricted software installation.

A key guiding principle of qLOOK is to keep metadata minimal but sufficient for complete reanalysis of experiments. This approach aligns with MIQE guidelines and FAIR data practices,

ensuring that stored qPCR data are transparent, reusable, and reproducible. By balancing simplicity with rigorous documentation, qLOOK provides a standardized framework suitable for individual laboratories, collaborative projects, and broader academic sharing.

Implementation and workflow

qLOOK is implemented as a modular, script-based workflow in R, consisting of three independent scripts—qLOOK_Module1_v1.0.R, qLOOK_Module2_v1.0.R, and qLOOK_Module3_v1.0.R. The modular architecture was intentionally chosen to allow users to execute only the relevant parts of the pipeline, such as re-running normalization or statistical analyses without repeating data extraction. Each module operates through graphical pop-up windows, requiring no coding skills or prior R experience (**Figure 1**).

qLOOK_Module1_v1.0.R performs data extraction and formatting. Upon selecting a folder containing .EDS files, the script automatically identifies and processes all available files— independently of cyclers model—and compiles the results into a structured spreadsheet, qLOOK_Data.xlsx, containing formatted Cq values and metadata. It also generates a reference gene stability report (qLOOK_RefGenes.xlsx) using ΔCq (4), GeNorm (5), and NormFinder (6) algorithms, and a log file (qLOOK_Summary.txt) that records all processing steps. Importantly, qLOOK_Module1_v1.0.R's data-parsing logic is modular by design. This allows the program to automatically recognize and process EDS archives with different internal structures. These include those generated by older software versions, such as 7500, 7500 Fast, StepOnePlus, Vii7, and QuantStudio6, as well as the newer Thermo Fisher DA2 real-time PCR software for QuantStudio3 and QuantStudio5 cyclers.

qLOOK_Module2_v1.0.R performs data normalization. After the user selects the qLOOK_Data.xlsx file, the script prompts for reference genes and a calibrator sample. The resulting files include qLOOK_Norm.xlsx (normalized and expression data) and qLOOK_Express.xlsx, which provides three linked sheets containing normalized values,

relative expression in linear format ($\Delta\Delta Cq$), and using the Livak method ($2^{-\Delta\Delta Cq}$; (7)). This module also generates distribution plots of $2^{-\Delta\Delta Cq}$ values and records all analytical steps, reference genes, and calibrator information in `qLOOK_Summary.txt`.

`qLOOK_Module3_v1.0.R` handles statistical analysis and data formatting for downstream applications. It imports `qLOOK_Express.xlsx`, requests a p-value threshold (default = 0.05), and performs one-way ANOVA and pairwise t-tests analyses. The results are stored in `qLOOK_ANOVA.xlsx` and `qLOOK_TTEST.xlsx`, including both unfiltered and significance-filtered comparisons. Additional outputs include a GraphPad-compatible file (`qLOOK_GraphPad.xlsx`) and a folder of automatically generated box plots (`qLOOK_Plots`). All computational steps are appended to the `qLOOK_Summary.txt` log, ensuring complete reproducibility.

Data structure and metadata model

The `qLOOK` system is designed around a minimal yet sufficient data model, ensuring that the essential elements required for qPCR data reanalysis are preserved in a standardized and transparent format. The output of the workflow centers on two primary files: `qLOOK_Data.xlsx` and `qLOOK_Summary.txt`.

`qLOOK_Data.xlsx` (**Table 1**) represents the structured data table generated from processed `.EDS` files. The first column lists sample identifiers, while each subsequent column corresponds to a target gene, with numerical entries representing `Cq` values. This matrix-style format allows straightforward visualization, reanalysis, and compatibility with standard statistical or plotting tools. Because the script can process large batches of `.EDS` files, users are advised to maintain consistent sample naming conventions and avoid spaces in names to ensure error-free handling in R. The table does not currently include experimental metadata (e.g., operator, instrument, or plate ID), as these parameters are already defined within MIQE

guidelines and are intended to be stored separately in laboratory documentation. However, the structure is readily extensible to accommodate future metadata layers if required.

qLOOK_Summary.txt (Table 2) functions as a comprehensive analytical log that records all key parameters influencing data interpretation. It includes script versions, timestamps, lists of processed .EDS files, instrument type, reference gene(s), calibrator sample(s), selected statistical thresholds, list of input files, and the names of all generated output files. Each computational step is chronologically documented in plain text, providing full transparency and reproducibility of the analytical workflow.

Together, qLOOK_Data.xlsx and qLOOK_Summary.txt constitute the minimal dataset required for complete reanalysis of qPCR experiments, in line with FAIR data principles. This format complements, rather than replaces, MIQE standards—translating them into a practical digital framework for ELN integration and long-term data stewardship

Reproducibility and interoperability

Reproducibility was a primary consideration in the design of qLOOK. Each analytical step is automatically recorded in qLOOK_Summary.txt, providing a transparent record of data extraction, normalization, and statistical analysis. The log includes all relevant parameters, such as script versions, reference genes, calibrators, and statistical thresholds, allowing users to reproduce identical results or apply alternative analytical settings using the same input data. Because all operations are performed through scripted procedures rather than manual data handling, qLOOK minimizes user-dependent variation and ensures consistency across analyses.

The separation between raw instrument output (.EDS files) and the standardized qLOOK_Data.xlsx format further supports reproducibility and long-term accessibility. Once extracted, data can be reanalyzed independently of the original proprietary software, reducing dependency on specific software versions or operating systems. By relying solely on open-

source R functions and generating plain-text or spreadsheet outputs, qLOOK maintains full transparency and platform independence.

Regarding interoperability, qLOOK was designed to align with FAIR and MIQE principles. The compact metadata structure allows straightforward incorporation into ELNs and compatibility with institutional or public data repositories. The modular code base also permits the inclusion of additional qPCR file formats or metadata fields without altering the core structure. Together, these features make qLOOK a practical and extensible tool for standardized digital storage and reproducible reanalysis of qPCR experiments.

Discussion and Conclusion

qLOOK provides a simple and standardized approach for storing, processing, and reanalyzing qPCR data independently of proprietary software. By retaining qLOOK_Data.xlsx and qLOOK_Summary.txt, users preserve the minimal information required to reproduce or extend analyses without reprocessing the original .EDS files, ensuring transparent data provenance and enabling future reanalysis. The current implementation requires R, RStudio, and RTools, but all necessary libraries are automatically installed and loaded, facilitating use by researchers without prior coding experience. A self-contained, R-independent Windows executable is under development to increase accessibility in environments with restricted administrative privileges. Although the current version focuses on Thermo Fisher/ABI .EDS files, the modular design of the scripts allows straightforward inclusion of additional file formats from other qPCR cyclers in the future. Overall, qLOOK represents a practical tool for minimal, standardized digital documentation of qPCR experiments, supporting reproducibility, transparency, and long-term data stewardship in molecular biology and biomedical research.

qLOOK Accessibility

The qLOOK scripts, ReadMe file, and license information are available for download from M.C.'s GitHub repository: <https://github.com/mircocastoldi>

Abbreviations

qLOOK; qPCR-LOG-boOK, MIQE; Minimum Information for Publication of Quantitative Real-Time PCR Experiments. qPCR; Quantitative real-time PCR. FAIR; Findable, Accessible, Interoperable, and Reusable, ELN; electronic laboratory notebooks.

References

1. *Bustin SA, Ruijter JM, van den Hoff MJB, Kubista M, Pfaffl MW, Shipley GL, et al. MIQE 2.0: Revision of the Minimum Information for Publication of Quantitative Real-Time PCR Experiments Guidelines. Clin Chem. 2025.*
2. *Bustin SA, Benes V, Garson JA, Hellemans J, Huggett J, Kubista M, et al. The MIQE guidelines: minimum information for publication of quantitative real-time PCR experiments. Clin Chem. 2009;55(4):611-22.*
3. *Wilkinson MD, Dumontier M, Aalbersberg IJ, Appleton G, Axton M, Baak A, et al. The FAIR Guiding Principles for scientific data management and stewardship. Sci Data. 2016;3:160018.*
4. *Silver N, Best S, Jiang J, Thein SL. Selection of housekeeping genes for gene expression studies in human reticulocytes using real-time PCR. BMC Mol Biol. 2006;7:33.*
5. *Vandesompele J, De Preter K, Pattyn F, Poppe B, Van Roy N, De Paepe A, et al. Accurate normalization of real-time quantitative RT-PCR data by geometric averaging of multiple internal control genes. Genome Biol. 2002;3(7):RESEARCH0034.*
6. *Andersen CL, Jensen JL, Orntoft TF. Normalization of real-time quantitative reverse transcription-PCR data: a model-based variance estimation approach to identify genes suited for normalization, applied to bladder and colon cancer data sets. Cancer Res. 2004;64(15):5245-50.*
7. *Livak KJ, Schmittgen TD. Analysis of relative gene expression data using real-time quantitative PCR and the 2^{(-Delta Delta C(T))} Method. Methods. 2001;25(4):402-8.*

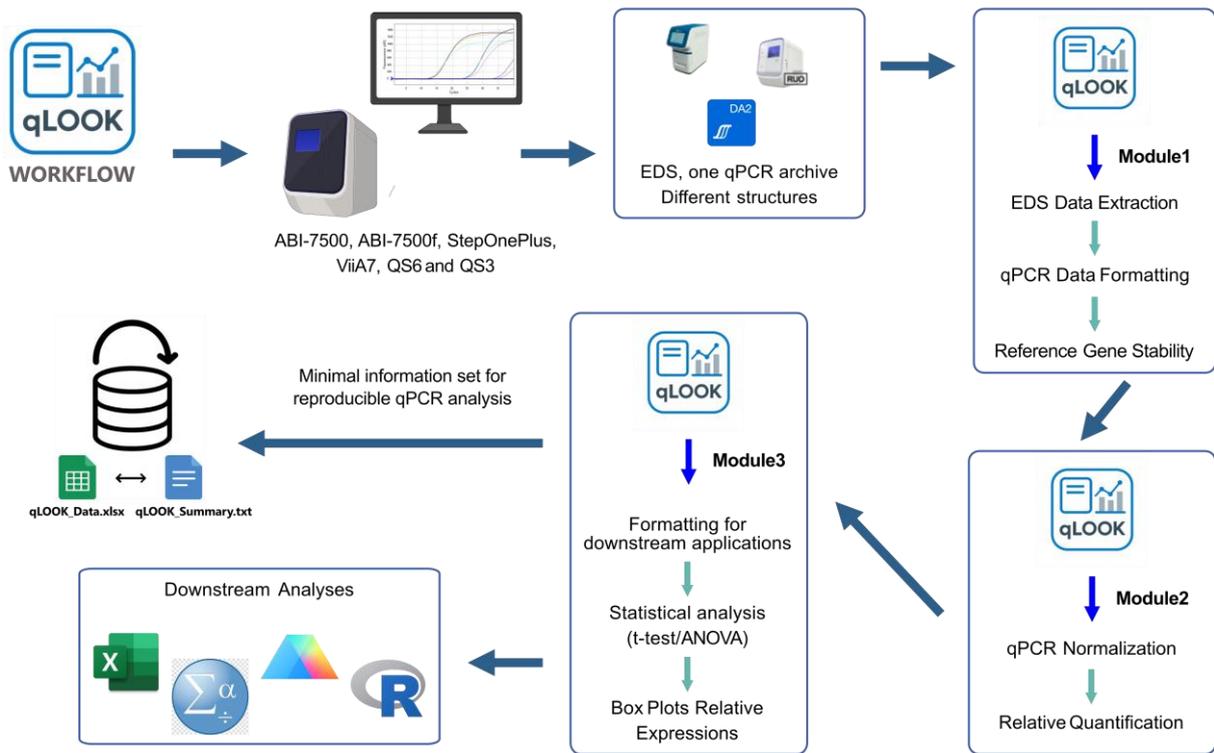


Figure 1 shows the schematic representation of the qLOOK Workflow.

Sample_Name	ANLN	AURKb	GAPDH	KIF11
phx_153_4h	22,3512	29,1157	20,6971	20,5673
phx_162_8h	21,0680	21,9792	21,4744	20,5262
phx_163_4h	22,8379	25,1766	20,6408	20,3407
phx_170_8h	21,8604	24,5755	21,1083	19,3488
phx_174_48h	14,9019	17,5523	19,8234	16,1102
phx_176_48h	14,6684	17,5083	18,9228	15,7772
phx_190_8h	21,2842	23,7515	19,9763	19,7595
phx_194_48h	14,4002	16,8717	19,1346	15,4296
phx_195_48h	15,2774	18,0843	20,1057	16,5263
phx_22_48h	15,9776	17,8521	20,8105	16,5250
phx_26_24h	16,5454	18,4588	20,7834	18,3250
phx_28_24h	18,4952	20,7846	20,6637	19,6882
phx_31_24h	17,5852	19,7550	19,4581	18,2496
phx_52_24h	21,3693	22,8459	20,3576	19,4385
phx_63_8h	19,9575	23,0572	21,2054	19,2469
phx_65_4h	21,4997	25,9957	21,9071	19,2268
phx_70_8h	25,2151	23,0949	20,5137	19,6649
phx_71_8h	20,7688	22,5058	20,2955	19,4323
phx_72_4h	21,8862	23,8683	21,4390	21,7104
phx_73_4h	23,0833	25,6393	21,7183	19,6814
sham_160_8h	34,2735	25,6100	21,2632	20,1316
sham_161_4h	21,6394	25,0499	20,5429	19,0115
sham_169_24h	22,9606	23,4582	20,2661	19,7673
sham_171_24h	20,8130	23,6791	20,7537	19,8000
sham_172_4h	21,1039	23,2030	20,1256	18,3395
sham_178_8h	23,6461	24,3788	21,8180	20,5194
sham_180_4h	20,9927	24,2251	20,9408	19,1501
sham_185_48h	20,9911	23,1441	21,6079	19,1118
sham_64_8h	21,6082	23,5504	20,3015	19,7384
sham_66_48h	20,3126	21,9276	21,2515	18,6223
sham_67_24h	23,6023	24,3852	20,2212	19,7568
sham_78_48h	20,4342	21,9305	20,8713	18,1395
sham_79_48h	20,5357	22,9319	20,1377	18,8086
sham_82_8h	22,1111	26,3318	20,8869	20,4121
sham_83_8h	22,1228	23,6353	20,2361	20,6040
sham_84_4h	21,4856	30,6392	21,3718	19,9124
sham_86_48h	19,4260	21,5193	20,2765	17,8837
sham_89_24h	22,1744	23,5963	19,8232	19,8377
sham_90_24h	19,8258	21,8817	20,3324	19,0025

Table 1 shows the content of the *qLOOK_Data.xlsx* file after being generated by the *qLOOK_Module1* script.

```

1  ----- qLOOK - qPCR-Log-boOK -----
2  ----- EDS Processing Log -----
3
4  --- 2025-10-15 13:42:21 ---
5  Script: qLOOK_Module1_v1.0.R (v1.0)
6  Script: qLOOK_RefGene v1.R (v1)
7  Original .EDS files processed:
8  - ANLN 05092025.eds processed as TXT (analysis_result.txt)
9  - AURKb 08092025.eds processed as TXT (analysis_result.txt)
10 - GAPDH 2025-09-10.eds processed as TXT (analysis_result.txt)
11 - KIF11 03092025.eds processed as TXT (analysis_result.txt)
12
13 If .TXT -> EDS was generated by ViiA7, StepOneplus, QuantStudio6, ABI7500, or ABI7500fast
14 If .JSON -> EDS was generated by QuantStudio3
15
16 Created files:
17 - qLOOK_qBASE_20251015_134207.xlsx
18 - qLOOK_qBASE_20251015_134207.txt
19 - qLOOK_Data_20251015_134207.xlsx
20 - qLOOK_RefGene_20251015_134207.xlsx
21
22
23 ----- qPCR Normalization and Relative Expression Log -----
24
25 --- 2025-10-15 13:44:00 ---
26 Script: qLOOK_Module2_v1.0.R (v1.0 (14-10-2025))
27 Input File: qLOOK_Data_20251015_134207.xlsx
28 Reference Gene(s): GAPDH
29 Calibrator Sample(s): sham_161_4h, sham_180_4h, sham_84_4h
30 Created Files:
31 - qLOOK_Norm_20251015_134236.xlsx
32 - qLOOK_Express_20251015_134236.xlsx
33 - Boxplots Folder: qLOOK_Norm_Plots_20251015_134353
34
35
36 ----- Statistical Analysis Log -----
37
38 --- 2025-10-15 13:46:14 ---
39 Script: qLOOK_Module3_v1.0.R (v1.0 (14-10-2025))
40 Input Excel: qLOOK_Express_20251015_134236.xlsx
41 p-value cutoff: 0.05
42 Created files:
43 - qLOOK_GraphPad_20251015_134607.xlsx
44 - qLOOK_ANOVA_20251015_134607.xlsx
45 - qLOOK_TTEST_20251015_134607.xlsx
46 Plots folder: qLOOK_Plots_20251015_134607
47 Plots:
48 - ANLN_boxplot.png
49 - AURKb_boxplot.png
50 - GAPDH_boxplot.png
51 - KIF11_boxplot.png
52 - Heatmap_20251015_134607.png
53
54

```

Table 2 shows the structure and content of the `qLOOK_Summary.txt` log file after the three `qLOOK` modules were run to analyze an experiment with four EDS files. In addition to logging the date, script version, and names of the input and output files, the `qLOOK` scripts log the metadata necessary for reanalyzing the experiment. This metadata includes the selected reference gene (`GAPDH`), the calibrator samples (`Sham_161_4h`, `Sham_180_4h`, and `Sham_84_4h`), and the *p*-value cutoff used to filter the results of the pairwise *t*-test (0.05).