



LabOS: The AI-XR Co-Scientist That Sees and Works With Humans

Le Cong^{1,2,*,**}, Zaixi Zhang^{3,*}, Xiaotong Wang^{1,2,*}, Yin Di^{1,2,*}, Ruofan Jin³, Michal Gerasimiuk^{1,2}, Yinkai Wang^{1,2}, Ravi K. Dinesh^{1,2}, David Smerkous⁴, Alex Smerkous⁵, Xuekun Wu^{2,6}, Shilong Liu³, Peishan Li^{1,2}, Yi Zhu^{1,2}, Simran Serrao^{1,2}, Ning Zhao^{1,2}, Imran A. Mohammad^{2,7}, John B. Sunwoo^{2,7}, Joseph C. Wu^{2,6}, Mengdi Wang^{3,**}

Affiliations:

¹ Department of Pathology, Department of Genetics, Stanford University School of Medicine, Stanford, CA, USA

² Institute for Stem Cell Biology and Regenerative Medicine, Stanford Cancer Institute, Stanford University School of Medicine, Stanford, CA, USA

³ Princeton AI Lab, Department of Electrical & Computer Engineering, Princeton University, Princeton, NJ, USA

⁴ School of Electrical Engineering and Computer Science, The Ohio State University, Columbus, OH, USA

⁵ Department of Bioengineering, University of Washington, Seattle, WA, USA

⁶ Division of Cardiology, Department of Medicine, Stanford University School of Medicine, CA.

⁷ Department of Otolaryngology-Head and Neck Surgery, Stanford University School of Medicine, Stanford, CA, United States

* Co-first and core contributing authors

** Corresponding authors: Le Cong (congle@stanford.edu), Mengdi Wang (mengdiw@princeton.edu)

Abstract: Modern science advances fastest when thought meets action. LabOS represents the first AI co-scientist that unites computational reasoning with physical experimentation through multimodal perception, self-evolving agents, and XR-enabled human-AI collaboration. By connecting multi-model AI agents, smart glasses, and human-AI collaboration, LabOS allows AI to see what scientists see, understand experimental context, and assist in real-time execution. Across applications—from cancer immunotherapy target discovery to stem-cell engineering—LabOS shows that AI can move beyond computational design to *participation*, turning the laboratory into an intelligent, collaborative environment where human and machine discovery evolve together.

Introduction

Science advances on two coupled fronts: computation that proposes and predicts, and experimentation that validates and reveals. Recent AI has transformed the computational front [1,2,3,4]—accelerating simulation, prediction, and design—yet the physical laboratory remains the bottleneck, where perception, coordination, and reproducibility still limit progress, and outcomes hinge on hard-to-transfer or tough-to-reproduce skills. Meanwhile, today’s “agentic AI” largely operates in the digital realm—planning experiments and synthesizing tools from text, data, and simulations without perceiving or acting in the dynamic lab. In parallel, robotic automation in laboratories can be powerful but mostly rule-based and bespoke—costly to retarget, challenging to relocate, and brittle to real-world variability.

LabOS addresses this gap through a unified human–AI collaborative intelligence platform that makes laboratories AI-perceivable and AI-operable. It integrates *agentic AI systems* for dry-lab

reasoning with extended reality(XR)-enabled, *multimodal interfaces* for human-in-the-loop wet-lab execution, creating an end-to-end framework that links hypothesis generation, experimental design, physical validation, and automated documentation.

The LabOS co-scientist adopts a multi-agent reasoning architecture—comprising planning, development, critique, and tool-creation agents—that collectively perform hypothesis generation, experimental design, data analysis, and adaptive improvement. The AI co-scientist is self-improving, continuously expanding its analytical capabilities via a “Tool Ocean” of modules autonomously generated from web search, scientific literature, and data. This self-evolving capacity enables the AI to solve novel research tasks via inference-time scaling.

In this paper, we focus on an end-to-end instantiation of LabOS for the biomedical domain and demonstrate state-of-the-art performance on leading biomedical reasoning benchmarks—including Humanity’s Last Exam (HLE): Biomedicine [5], LAB-Bench: DBQA, and LAB-Bench: LitQA [6]—while closing the loop from dry-lab planning to wet-lab execution.

On the physical side, LabOS connects AI reasoning directly to the laboratory via AR/XR smart glass interfaces and real-time multimodal sensing. Researchers wearing XR glasses receive adaptive, context-aware guidance from the AI agent—step-by-step instructions, error detection and correction cues, and gesture or voice interactions for sterile workflows.

To enable the AI co-scientist to “see” in the lab, we collected >200 egocentric video sessions from researcher-worn cameras/glasses during real experiments, assembling these into the LabSuperVision (LSV) benchmark for evaluating AI models’ lab perception and reasoning capabilities. Because leading AI models struggled on this benchmark, we post-trained a lab-specialized Vision-Language-Model (VLM) using a combination of public experimental videos, in-house recordings, and expert annotations. The resulting model, namely the LabOS VLM, is able to decode visual input from XR glasses and align the visual embedding with a language model to interpret and reason about lab video scenes. The model demonstrates markedly improved visual reasoning capability in scientific settings, enabling LabOS to monitor actions, detect deviations, verify results, and synchronize multimodal data streams with reference protocols—allowing the AI to perceive, understand, and act/co-pilot within real laboratory as an AI co-scientist.

LabOS also supports 3D/4D spatial modeling of laboratory workflows. These digital twins capture spatial and temporal relationships between instruments, samples, and human actions, enabling replay, “what-if” analysis, and simulation-based training. The resulting spatial grounding provides the foundation for safe, reproducible, and transferable laboratory automation.

For real world validation, we demonstrate LabOS’s ability in three biomedical research studies: cancer immunology, mechanistic research, and stem cell engineering (**Fig. 1c**). In the first study, we sought to uncover cancer immunotherapy targets that could boost natural killer (NK) cell killing of tumors. We used LabOS to generate hypotheses, perform target identification via multi-step reasoning and analysis, and the AI agent nominated CEACAM6 as a putative target, which were validated in wet-lab NK-tumor killing assay. In the second study, we used LabOS to perform mechanistic research, showcasing identification of a gene regulating cell fusion, ITSN1. In the third study, two researchers wore smart glasses during stem-cell engineering projects and interacted with the AI co-scientist. The latest LabOS streams egocentric data to the server and invokes the VLM agent every ~5-10s. In a gene-knockout experiment, LabOS monitored the

workflow, provided step-level guidance, and could flag operation deviations—such as when the researcher did not follow sterile techniques, or, wrong reagent incubation time versus gold-standards. In a lentiviral transduction in human iPSCs experiment, LabOS was used to automatically record and digitalize workflows by an experienced researcher, which can be used for coaching a novice researcher to complete the experiment, demonstrating its potential in capturing and transferring advanced human skills.

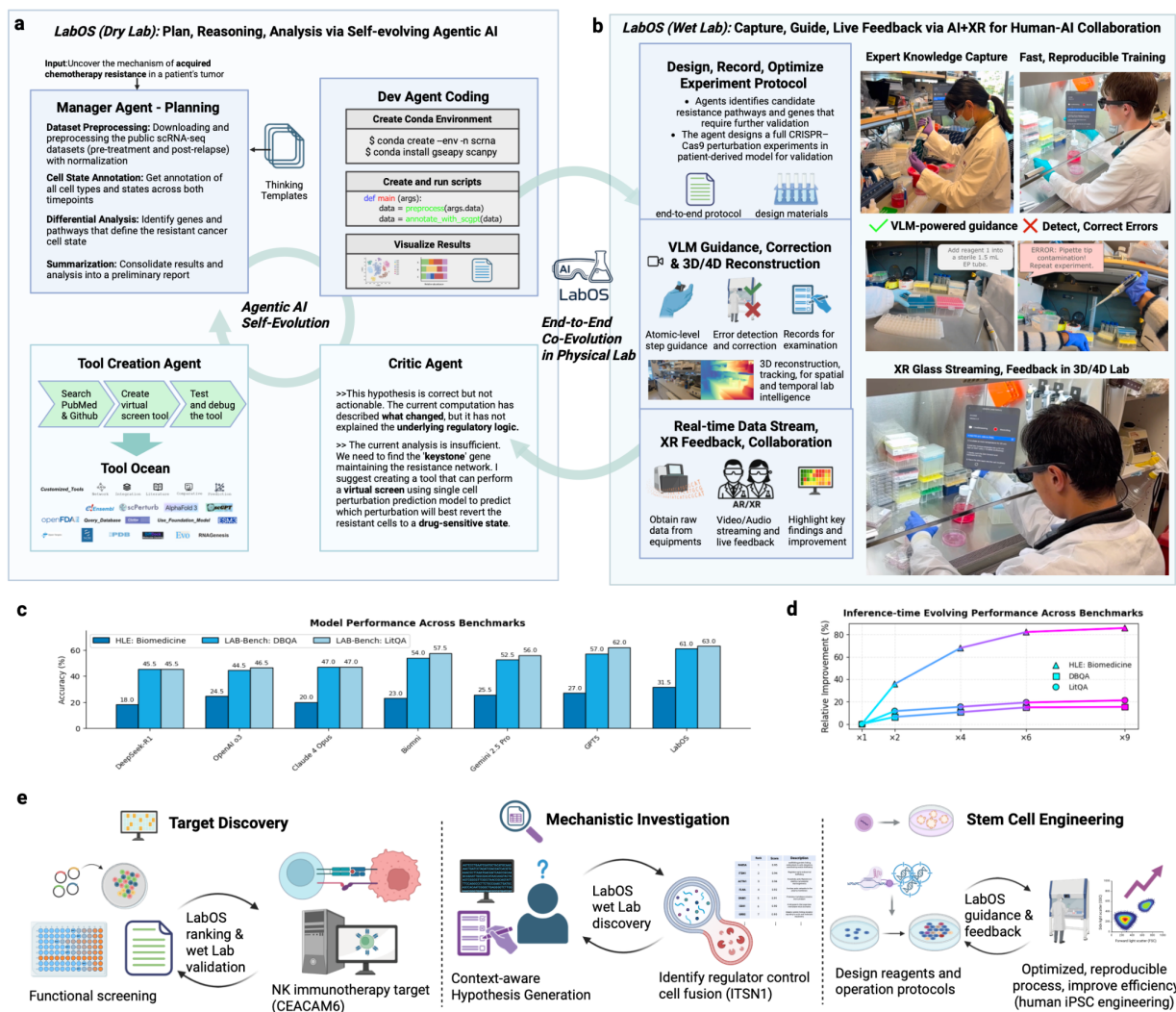


Figure 1. LabOS: Multi-Modal Human-AI Collaboration in Science Laboratory. LabOS comprises a self-evolving agentic AI for dry lab tasks and an XR interface for human-in-the-loop wet lab execution, creating an end-to-end system for lab research. **a, Dry Lab:** Agentic AI as computational core of LabOS. The system employs a multi-agent framework, including a Manager Agent for planning, a Dev Agent for coding and execution, and a Critic Agent for evaluation and refinement. The agents continuously improve analytical capabilities. A Tool Creation Agent expands the system's functions by generating new tools from sources like PubMed, which are then stored in a shared Tool Ocean. This component automates complex data analysis. **b, Wet Lab:** Capture, guide, live feedback via AR/XR for human-AI collaboration in the physical laboratory. A scientist wearing AR/XR glasses receives real-time guidance from the AI. A specially trained Vision-Language Model (VLM) monitors the procedure, providing on-the-fly error detection and correction to ensure correctness and reproducibility. **c, Results on leading benchmarks—HLE: Biomedicine, LAB-Bench: DBQA, and LAB-Bench: LitQA—**show LabOS outperforming frontier LLMs/agents in biomedical reasoning tasks. **d, Self-evolving agent's performance scales with inference-time compute.** **e, Use Cases.** The first case is drug target discovery: the agentic AI analyzes functional screening data to identify and rank NK cancer immunotherapy targets. Secondly, for mechanistic investigation, LabOS generates testable hypotheses that are then validated by a human scientist to identify a cell fusion regulator. Thirdly, LabOS enabled copiloted, reproducible processes for complex experiments like stem cell engineering. **All photographs shown are of the authors.**

Together, these features make LabOS a true AI co-scientist—the first multimodal human–AI system to unify dry-lab reasoning and wet-lab execution in a single, adaptive framework. By giving AI the ability to think with us and work alongside human scientists—seeing what we see, checking what we do, and learning from every run—LabOS turns the lab into a dynamic feedback loop. It moves us toward autonomous, self-improving discovery, where human intuition and machine rigor co-evolve to accelerate breakthroughs.

Results

1. LabOS Overview

LabOS links agentic AI for dry-lab reasoning with XR-enabled, human-in-the-loop wet-lab execution into an end-to-end workflow from design to validation. A multi-agent core—Manager (planning), Developer (coding/execution), and Critic (evaluation)—plus a Tool-Creation module that auto-extends a shared Tool Ocean nominates targets, plans experiments, runs analyses, and continually improves [7,8] (**Fig. 1a**).

For multimodal human-AI collaboration in the wet-lab module, protocols can be generated by or input into LabOS, then streamed to the connected XR glasses. The interface on XR glasses (i) renders stepwise protocol in an Unity/Android application, (ii) verifies physical actions from the first-person video stream by invoking an embedded VLM for visual reasoning, and (iii) returns context-aware feedback in real time (**Fig. 1b**). All streams are time-stamped and logged with metadata for automated documentation.

2. LabOS Supports Self-Evolving AI Agent for Biomedical Reasoning

The dry-lab module of LabOS builds on and expands the STELLA self-evolving agent framework [8], a multi-agent reasoning system for biomedical research that integrates planning, development, critique, and tool creation. **Fig. 1(a)** illustrates the agentic architecture. The Manager/Planner Agent decomposes scientific objectives into structured modules—candidate molecules, reagents, procedures, materials lists, instrument settings, and quality control checkpoints. The Developer Agent executes these steps by generating and running Python code for complex bioinformatics analyses, while the Critic Agent evaluates intermediate results and refines the workflow, forming an iterative reasoning loop.

This AI co-scientist learns and improves from every problem it solves, continuously enhancing its technical abilities by proactively expanding reasoning strategies and creating new tools [7,8]. Two mechanisms drive its continuous self-improvement. First, a Template Library of successful reasoning workflows is dynamically updated, allowing the system to generalize from prior solutions. Second, a Tool Ocean maintains and expands a repository of analytical tools/codes, databases, and APIs. The Tool-Creation Agent autonomously identifies, tests, and integrates new resources as needed. This architecture enables LabOS to plan and solve complex research tasks efficiently while continually enhancing its reasoning and analytical capabilities.

To evaluate the efficacy of LabOS research agent, we test-ran it on a suite of challenging biomedical reasoning benchmarks. Our results show that LabOS consistently establishes a new state of the art, achieving top accuracy scores of approximately 32% on *Humanity's Last Exam: Biomedicine*, 61% on *LAB-Bench: DBQA*, and 65% on *LAB-Bench: LitQA*, outperforming the

next-best models by up to 8% (**Fig. 1c**). Crucially, its performance systematically improves with continued use and test-time scaling, providing direct evidence of its self-evolving design (**Fig. 1d**). The self-improving framework enables LabOS to learn and grow like a human scientist, dynamically scaling to meet the ever-expanding challenges of biomedical discovery.

3. Training LabOS AI Co-Scientist To See and Reason in Physical Labs

To enable LabOS to see, understand, and reason in a physical lab environment, we sought to post-train a Vision-Language-Model (VLM) on a broad set of lab research videos. A VLM is a multimodal AI model that jointly learns from visual and textual inputs, allowing it to connect what it *sees* with what it *reads* or *describes*. It typically combines a vision encoder that processes images or video frames with a language model that interprets and generates text. By aligning these two modalities in a shared representation space, the VLM can recognize lab scenes, interpret actions, and reason about experimental workflows and outcomes.

3.1 LabSuperVision (LSV): Benchmarking AI for Scientific Visual Reasoning

We asked researchers to wear XR glasses or action cameras when they operate in their laboratory research to collect real-world data. Then we processed and annotated the collected data for evaluating AI models in lab research.

We assembled **LabSuperVision (LSV)**, an expert-annotated laboratory video dataset, designed for lab operation video understanding and reasoning. LSV comprises >200 high-quality video sessions (typically 2–10 min; up to 45 min) recorded by 7 researchers across diverse lab scenes—bench, tissue culture room, and instrument bays—capturing realistic operations and movement between spaces. Each session is associated with an expert-generated gold-standard protocol. Then, human scientists performed corresponding experiments wearing cameras or smart glasses, recording all details. For annotation, a team of five human experts annotate each video session with: (1) Step segments with start/stop times aligned to the gold-standard protocol. (2) Error and issue events labeled by type (e.g., sterile breach, step mismatch, timing deviation). (3) Critical parameters, materials and reagents when applicable (**Fig. 3**).

With the LSV benchmark, we evaluate whether leading commercial AI models can interpret and troubleshoot scientific procedures from lab experiment videos. We tested four leading AI models spanning both open and proprietary foundation models: Gemini 2.5 Pro, Cosmos-Reason-1, Qwen 2.5-VL-7B, and GPT-4o (**Fig. 3c**). When prompted with a new video, each model is tasked with one of the two: (1) Protocol alignment, where the model needs to generate a stepwise protocol describing procedural actions and parameters, (2) Issue identification, where the model needs to troubleshoot based on the video and identify any deviation from gold-standard or any handling errors if present. The model outputs are then compared against the ground truth, using two methods: a total of five human experts compared VLM outputs vs. gold-standard protocol and error labeling following a 0-5 scoring (5 being the highest for perfect alignment); in parallel, we used GPT-5 as a comparator to examine and score outputs using the same rubrics.



Figure 2. LabOS-VLM for visual reasoning in laboratories. **a**, LabSuperVision (LSV) benchmark dashboard with egocentric and fixed-view lab videos paired to gold-standard protocols, including notes, issues or errors labeled by human experts. **b**, LSV benchmarking pipeline in four stages: (1) Data collection—multi-modal recordings from diverse facilities via fixed cameras and XR smart glasses; (2) Data processing—expert curation with light automation, aligned to reference protocols and parameters; (3) Dataset assembly—compressed corpus of 200+ experiment-session videos; (4) Evaluation—human and GPT-5 assessments of model performance on LSV. **c**, Benchmarking AI models on LSV via human experts (n=5) and GPT-5. Scores range 0–5 (5 = perfect). **d**, LabOS-VLM post-training using FineBio, JoVE, and LSV datasets via supervised fine-tuning (SFT) and reinforcement learning (GRPO). **e**, LabOS-VLM family outperform baselines on protocol generation quality (left) and error-detection rate (right) on LSV. **f–g**, Real-world testing of LabOS-VLM: (f) error detection/correction in a tissue culture experiment; (g) context-aware AI generated step-by-step instructions grounded in visual understanding.

Our findings show that leading AI models struggle to understand fine-grained laboratory research workflows: the top-performing model, Gemini-2.5Pro, scored only 2.86 out of 5 in protocol alignment, moderately better than open-source NVIDIA Cosmos-1 which scored 2.24; for issue/error identification, leading models like Gemini, GPT4o only managed to score ~2 out of 5 (**Fig. 2c**). This reflects the gap between open-world foundation models' capabilities and what is needed for scientific research in physical labs. Our results surface where today's AI models have some success (protocol alignment) and where they have substantial limitations (error detection/correction).

3.2 LabOS-VLM: Vision-Language-Model Trained for Visual Reasoning in Laboratories.

Given the aforementioned gap, we sought to train a vision–language model for visual reasoning in lab environments. For this purpose, we assemble an experimental video collection, including FineBio (expert-annotated wet-lab videos), JoVE (standardized procedure videos), and LSV. We split the datasets 80/10/10 for training, validation, and held-out testing.

Using Qwen-VL as the base model, we performed post-training by supervised fine-tuning (SFT) with LoRA on paired video–text examples and then reinforcement finetuning to improve visual reasoning. In the latter reinforcement learning step, we used Group Relative Policy Optimization (GRPO) [9] with LoRA updates, where the model rolled out a group of multiple candidate responses per prompted scenario and received rewards designed to evaluate the candidate outputs. The reward is designed to be rule-based to account for lab procedural accuracy, safety compliance, and experimental detailed level; as well as to account for relative rewards within each group to favor expert-consistent reasoning. By using this SFT->RL pipeline, we adapted the base model and obtained the LabOS-VLM (7B, 32B, 72B, 235B; **Fig. 2e–f**) family.

Across model scales, the LabOS-VLM family consistently outperforms the base model on lab video reasoning tasks (measured via protocol generation quality and error detection accuracy). On the held-out subset of evaluation data, LabOS VLM-235B achieves >90% accuracy in error detection accuracy, outperforming Claude Opus-4.1, GPT-5, and Gemini 2.5 Pro on all evaluated metrics, establishing LabOS VLM as a strong module for multimodal lab intelligence and AI-assisted experimentation (**Fig. 2e**). Our SFT->RL pipeline has proved to enhance the AI model's ability in step-wise visual reasoning, error detection and correction, which is a basis for building AI co-scientist to see and work with researchers real-time in scientific settings.

We next evaluated the fine-tuned LabOS-VLM on egocentric videos from real experiments. First, we supplied a series of videos about delivering CRISPR gene-editors to human cells via transfection. The model successfully distinguished correct vs. incorrect operations, pinpointing the specific errors that the human operator made across 2 distinct issues (**Fig. 2f**). Second, we provided experiment records, from a multi-step setup on preparing Cas9 RNA complex for gene-editing along with a reference protocol, to the LabOS co-scientist. The VLM-powered AI agent recognized each step, generated step-wise guidance, matched actions to the protocol, issued context-aware warnings, and suggested the next action. These results validate LabOS VLM capabilities in authentic wet-lab settings.



Figure 3. LabOS on XR glasses enables spatially-grounded human-AI collaboration in physical laboratories. **a**, Live streaming from XR glasses to the LabOS VLM server enables real-time scene understanding, feedback, and interactions between human and AI agents. **b**, Deployment of LabOS AI+XR system across lab settings. **c**, Live action feed from the wearer's perspective. **d**, LabOS AI+XR provides guidance and summary on the lab operation. **e**, Detected deviations trigger inline error prompts and suggested corrections, to mitigate human researcher's oversight. **f**, Feed-forward Gaussian splatting supports camera tracking and multi-view, metric-depth 3D reconstruction, enabling object localization and spatio-temporal reasoning for scientific workflows. **All photographs shown are of the authors.**

4. Extended-Reality (XR) Glasses Enable Real-Time Human-AI Collaboration

Human–AI interaction in LabOS is mediated through the XR interface. Researcher engages with the AI co-scientist primarily via the LabOS XR glasses, where the glasses live-stream what the human scientist sees and hears while reconstructing the 3D lab environment. Modern XR hardware supports interface rendering, gesture recognition, and running a Unity/Android application. Data from the glasses are streamed to a local GPU server (or the cloud) for real-time agentic inference. The server receives short video segments (e.g., 5–10s), forwards them to the LabOS AI agent for analysis and reasoning (input at 4 frames/s for videos), and returns a structured JSON output to the XR glasses. The XR app parses this JSON message and provides real-time visual and audio feedback to researchers at the bench (**Fig. 3a-e**). We tested both AR/XR glasses and VR/XR headsets, while both supported the above functions, we chose to launch initial deployment of LabOS using AR/XR glasses as they achieved key specifications needed for human-centered intelligent lab: open-peripheral, light-weight hardware (less than 3-ounces/85-grams) allows convenient wearing in labs where goggles are already used; 2+ hours battery life for extended operation (via a wireless neckband or mobile battery); display brightness over 1200+ Nits sufficient for in-door labs; 6DoF and hand gesture support for 3D-aware human-AI interactions.

LabOS also supports 3D modeling of lab environments and workflows (**Fig. 3f**), via videos captured on researcher-worn smart glasses — optionally augmented with multiview cameras — using state-of-the-art 3D/4D reconstruction algorithms. In our experiment, we utilize MapAnything [10] for multiview environment + egocentric for tracking, camera positioning, depth maps, and point cloud reconstructions, and incorporate 3D Gaussian splatting. Gaussian splatting models scenes as sets of millions of Gaussian distributions whose parameters are inferred from input 2D images captured in different camera positions and (and time if needed), enabling photorealistic, temporally consistent reconstruction. 4DLangSplat [11] can further help produce a time-aware, semantically indexable 3D environment, to support object-centric tracking.

Through remote supervision and guidance, the LabOS AI Co-Scientist boosts reproducibility by standardizing experiment capture and logging non-fatal deviations and context variables. Expert-run recordings from XR glasses also serve as instructional material and training data to improve LabOS agents and downstream models.


5. LabOS In Action: Human-AI Collaboration for Biomedical Discovery

5.1 LabOS for cancer immunotherapy target discovery


The challenge in cancer immunology is to identify regulators of immune evasion that mediate tumor resistance to cytotoxic lymphocytes such as natural killer (NK) cells. Traditional screening approaches are limited by throughput and rely on human expertise for downstream analysis. We prompted LabOS to investigate genes regulating sensitivity and resistance to NK cell-mediated killing of melanoma cancer cells based on a functional screen.

Using a functional CRISPR activation (CRISPRa) screen in A375 melanoma tumor cells, treated with or without primary human NK cells, LabOS AI co-scientist autonomously identified and refined candidate regulators of NK-mediated cytotoxic resistance [12]. Within a wet-lab-in-the-loop pipeline, the LabOS AI agent dynamically re-ranked genes through iterative functional enrichment analysis, revealing CEACAM6 as a top regulator of NK resistance.


a. LabOS accelerates functional screening for cancer immunotherapy target discovery

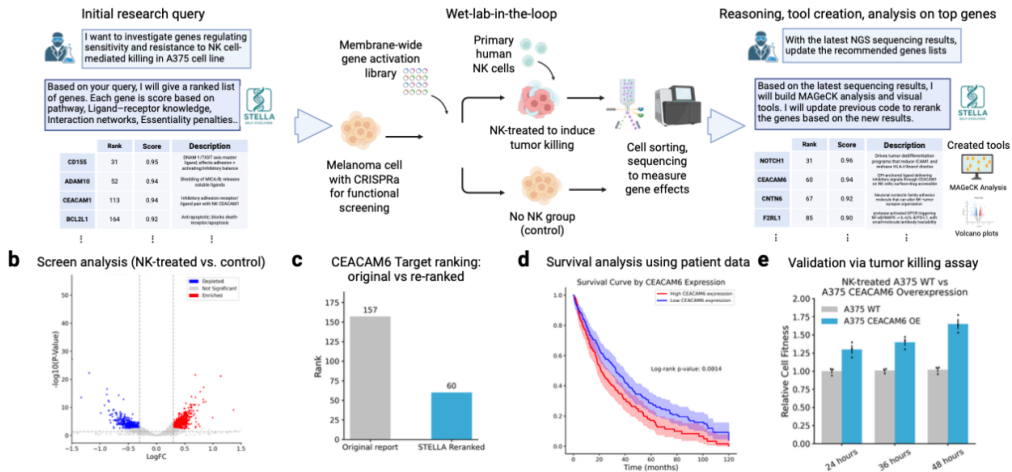

f. LabOS facilitates mechanistic investigation of genes that control cell fusion

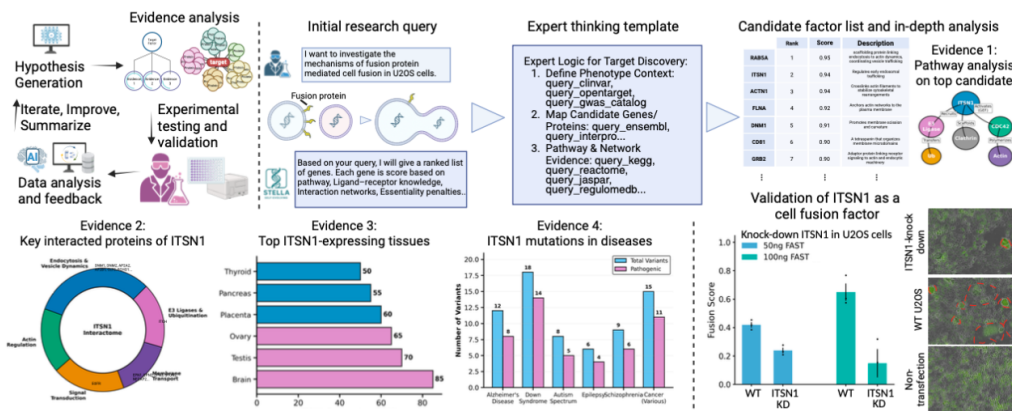

g. LabOS enables live guidance, error detection to safeguard multi-step experiments

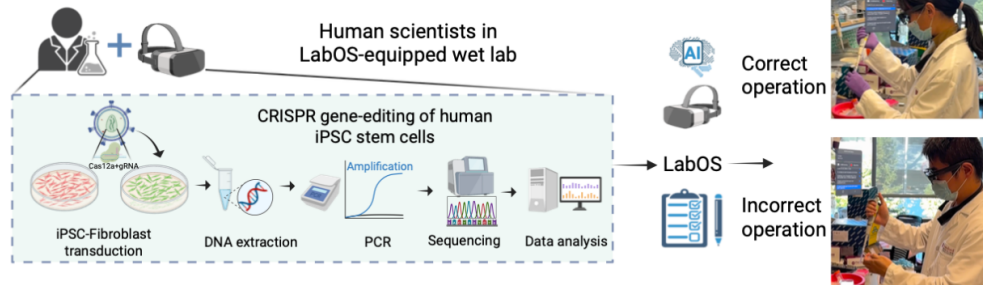

h. LabOS learns from experts to train junior scientist to perform stem cell experiments

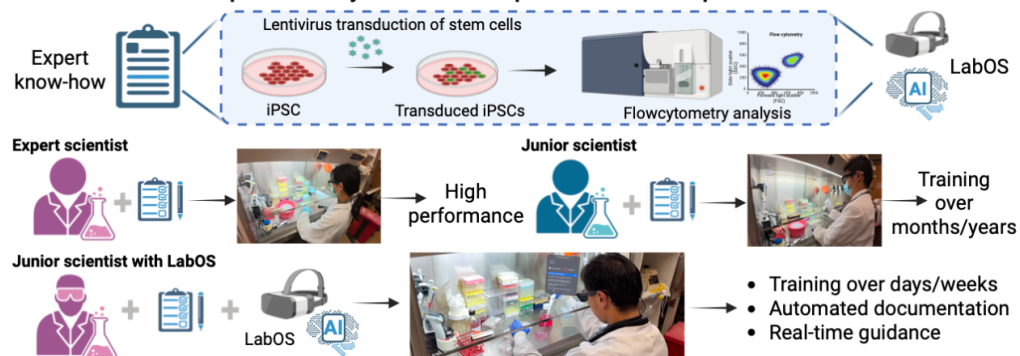


Figure 4. LabOS applications in target discovery, mechanistic investigation, and stem cell research. a-e, Functional screening study for Natural Killer (NK) immunotherapy target identification. LabOS was prompted with the research task to identify regulators of tumor resistance to NK killing, based on results from a CRISPR activation screen in A375 melanoma cells co-cultured with primary human NK cells. **b-e**, AI-generated screen analysis reveals differential gene enrichment upon NK treatment (b), and the AI agent re-ranked key target—*CEACAM6*—moving it from low to top-ranked

position (c) and automated survival analysis using cancer patient data, stratified by CEACAM6 expression (d). Wet-lab functional validation confirms enhanced NK resistance upon CRISPRa activation of CEACAM6 (e). **f**, LabOS applies to mechanistic investigation of cell fusion factors. LabOS proposes ITSN1 as a regulator of cell-cell fusion, providing its found evidence and reasoning trajectories. Experimental knock-out of top candidate gene ITSN1 confirmed this gene's impact on cell fusion in the wet lab, validating the AI co-scientist's hypothesis. **g**, LabOS provides live copilotting to researchers in stem cell gene-editing experiments. The AI-XR agent enables live guidance, error detection to safeguard multi-step experiments, where scientists can use LabOS to monitor and track issues automatically. **h**, The AI-XR agent can learn advanced skills in experiments and document them such as the lentiviral transduction of human iPSC stem cells. **All photographs shown are of the authors.**

AI reasoning further generated explainable evidence by automatically performing survival analyses on *The Cancer Genome Atlas* (TCGA) datasets, stratified by CEACAM6 expression, thereby linking functional screening results to patient outcomes. Experimental validation using individual CRISPRa perturbation confirmed that CEACAM6 activation significantly increased tumor resistance to NK killing (**Fig. 4a-g**).

5.2 LabOS guides mechanistic investigation of cell fusion regulator

Beyond functional screening, LabOS extends to mechanistic hypothesis generation and validation. In the next study, we prompt LabOS with the research task to identify key genes that control cell-cell fusion. Cell fusion is a fundamental biological process, key to basic biology (e.g. muscle development, viral infection all require cell fusion) and translational technology (e.g. cell fusion is a key step mediating efficient gene delivery).

LabOS proposed and ranked candidates using pathway enrichment, interaction priors, and functional evidence. The AI co-scientist prioritized *ITSN1* as a key regulator of cell-cell fusion, with automated evidence generation. Then human researchers used CRISPR interference (CRISPRi) coupled with cell fusion assay (induced by the fusogenic protein FAST [13]) in U2OS cells, to measure if *ITSN1* perturbation will impact cell fusion phenotype. Researchers further performed quantitative imaging and cell-based assay and observed significant inhibition of cell fusion upon *ITSN1* knockdown in the wet lab, confirming the role of *ITSN1* (**Fig. 4f**).

5.3 LabOS copilots researchers in complex stem-cell experiments

Reproducibility in advanced wet-lab domains is hindered by tacit expertise: critical steps are encoded in human memory and lab-specific “know-how,” not necessarily in written protocols. For stem cell engineering, as an example, small deviations in the timing, reagent handling, or cell density and cell state assessment can drive large outcome variation, impeding transfer of skills across operators and laboratories.

First, we used LabOS with XR glasses to guide researchers through CRISPR gene-editing of human stem cells for disease-modeling workflows [14]. Via XR glasses, LabOS agents captured all details as scientists perform gene knock-out experiments in human induced pluripotent stem cells (iPSCs)-derived cardiac fibroblast, useful for modeling diseases such as heart fibrosis [15]. LabOS copiloted the multi-step procedure with visual reasoning, precisely interpreting bench actions and flagging issues in real time (**Fig. 4g**).

Second, because efficient gene delivery underpins construction of engineered stem-cell lines for perturbation and drug screening, we used LabOS with XR glasses to guide and automatically document lentiviral transduction experiments in human iPSCs from expert-level scientists (**Fig. 4h**). LabOS then could monitor and help train junior scientists on this complex experiment.

Here, LabOS can act as an AI tutor: recording expert practice, digitizing key parameters, and coaching a novice to expert-level performance without requiring side-by-side training or extended period of trial-and-error. Through real-time, multimodal human-AI interactions, the LabOS system can provide context-aware guidance, validation, and performance tracking.

Summary

LabOS prototypes what an AI co-scientist can be: a system that sees, reasons, and helps run the lab. By pairing AI agents with real-time, XR-guided human–AI interaction and data-driven reasoning, it enables faster discovery, reproducible training, and precise operation. Across use cases—hypothesis generation, automated documentation, error correction, rapid skill transfer, iPSC experiment guidance, and insights into NK–tumor pathways and cell-fusion regulators—LabOS turns the lab into an adaptive, collaborative workspace where human scientists and AI work side-by-side to accelerate discovery, generate reproducible science, and improve together.

Acknowledgement

We extend our sincere gratitude to the NVIDIA, in particular the XR/AR-VR, VLM/Agentic-AI, edge computing, and business partnership team for their invaluable collaboration on the VLM fine-tuning development, agentic AI integration, and XR live streaming, as well as for GPU server support. We would like to acknowledge we procured the AR/XR glass hardware from Viture, and thank the Viture team for support on wireless streaming, Unity/Android app support. We would like to acknowledge we procured the VR/XR glass thanks to Dr. Ze Yuan from the UReality team and their software support. We also would like to thank the Nebius team for GPU compute and infrastructure support. The results on cancer patient survival analysis here are in whole or part based upon data generated by the TCGA Research Network: <https://www.cancer.gov/tcga>.

Author contributions

L.C. and M.W. conceived and supervised the project. L.C. led the research plan/organization, hardware/XR solution build, study design and paper writing. AI & Data: Z.Z., R.J. co-led the STELLA AI agents and benchmarking. X.W. led the data collection, annotation and human evaluation efforts, contributed to setting up XR glasses and live experimentation. Z.Z., Y.W., M.G. contributed to VLM training and benchmarking. D.S., S.L., A.S. co-led 3D/4D modeling, object recognition and tracking of lab workflows. D.Y., X.W., P.L., Y.Z., S.S., N.Z. contributed to data collection, scientific annotation, and human evaluation of VLMs. Experiment: X. W and R.D. co-led NK experiments, to which J.W. also contributed and supervised. D.Y. led the fusogen experiment and co-lead stem cell experiments with X.W., J.C.W. co-supervised the stem cell experiments.

Reference

1. Jumper, J. et al. Highly accurate protein structure prediction with AlphaFold. *Nature* 596, 583–589 (2021).
2. Wang, H. et al. Scientific discovery in the age of artificial intelligence. *Nature* 620, 47–60 (2023).
3. Baker, B. et al. Emergent autonomous scientific research capabilities of large language models. *Nat. Mach. Intell.* 6, 876–887 (2024).
4. Qu, Y. et al. CRISPR-GPT for agentic automation of gene-editing experiments. *Nat. Biomed. Eng.* (2025).
5. Phan, L. et al. Humanity's Last Exam (HLE). arXiv preprint arXiv:2501.14249 (2025).
6. Laurent, J. M., Janizek, J. D., Ruzo, M. et al. LAB-Bench: Measuring Capabilities of Language Models for Biology Research. arXiv preprint arXiv:2407.10362 (2024).
7. Qiu, J. et al. Alita: Generalist agent enabling scalable agentic reasoning with minimal predefinition and maximal self-evolution. arXiv preprint arXiv:2505.20286 (2025).
8. Jin, R. et al. STELLA: self-evolving LLM agent for biomedical research. arXiv preprint arXiv:2507.02004 (2025).
9. DeepSeek Team. Group Relative Policy Optimization (GRPO) for reinforcement learning of reasoning. arXiv preprint arXiv:2504.09374 (2025).
10. Keetha, N., Müller, N., Schönberger, J., Porzi, L., Zhang, Y., Fischer, T. & Kotschieder, P. MapAnything: Universal feed-forward metric 3D reconstruction. arXiv preprint arXiv:2509.13414 (2025).
11. Li, W., Zhou, R., Zhou, J., Song, Y., Herter, J., Qin, M. & Pfister, H. 4D-LangSplat: 4D language Gaussian splatting via multimodal large language models. In *Proc. CVPR*, 22001–22011 (2025).
12. Dinesh RK, Wang X, Mohammad IA, Gunasekaran P, Stikloraitis K, Villafuerte JR, Rao A, Hernandez-Lopez RA, Sunwoo JB, & Cong L. Surfaceome CRISPR activation screening uncovers ligands regulating tumor sensitivity to NK cell killing. *bioRxiv*. 2025:2025-09 (2025).
13. Duncan, R. Fusogenic reoviruses and their fusion-associated small transmembrane (FAST) proteins. *Annu. Rev. Virol.* 6, 341–363 (2019).
14. Takahashi K, Yamanaka S. Induction of pluripotent stem cells from mouse embryonic and adult fibroblast cultures by defined factors. *Cell*. 126(4):663-76 (2006).
15. Zhang H, Thai PN, Shivnaraine RV, Ren L, Wu X, Siepe DH, Liu Y, Tu C, Shin HS, Caudal A, Mukherjee S. Multiscale drug screening for cardiac fibrosis identifies MD2 as a therapeutic target. *Cell*. 187(25):7143-63 (2024).