

SignalMC-MED: A Multimodal Benchmark for Evaluating Biosignal Foundation Models on Single-Lead ECG and PPG

Fredrik K. Gustafsson^{1,2} Xiao Gu¹ Mattia Carletti¹ Patitapaban Palo¹
 David W. Eyre^{2,3,4} David A. Clifton^{1,2,5}

¹Department of Engineering Science, University of Oxford

²NIHR Health Protection Research Unit in Healthcare Associated Infections and Antimicrobial Resistance

³Big Data Institute, Nuffield Department of Population Health, University of Oxford

⁴NIHR Oxford Biomedical Research Centre

⁵Oxford Suzhou Centre for Advanced Research, University of Oxford

fredrik.gustafsson@eng.ox.ac.uk, xiao.gu@eng.ox.ac.uk, david.clifton@eng.ox.ac.uk

Abstract

Recent biosignal foundation models (FMs) have demonstrated promising performance across diverse clinical prediction tasks, yet systematic evaluation on long-duration multimodal data remains limited. We introduce SignalMC-MED, a benchmark for evaluating biosignal FMs on synchronized single-lead electrocardiogram (ECG) and photoplethysmogram (PPG) data. Derived from the MC-MED dataset, SignalMC-MED comprises 22,256 visits with 10-minute overlapping ECG and PPG signals, and includes 20 clinically relevant tasks spanning prediction of demographics, emergency department disposition, laboratory value regression, and detection of prior ICD-10 diagnoses. Using this benchmark, we perform a systematic evaluation of representative time-series and biosignal FMs across ECG-only, PPG-only, and ECG + PPG settings. We find that domain-specific biosignal FMs consistently outperform general time-series models, and that multimodal ECG + PPG fusion yields robust improvements over unimodal inputs. Moreover, using the full 10-minute signal consistently outperforms shorter segments, and larger model variants do not reliably outperform smaller ones. Hand-crafted ECG domain features provide a strong baseline and offer complementary value when combined with learned FM representations. Together, these results establish SignalMC-MED as a standardized benchmark and provide practical guidance for evaluating and deploying biosignal FMs.

Foundation models (FMs) for time-series and biosignals have recently emerged as a promising paradigm for scalable representation learning across diverse downstream tasks [7, 11, 13, 17, 29]. In cardiovascular medicine, large-scale pretraining on electrocardiogram (ECG) and photoplethys-

mography (PPG) signals offers the potential to learn physiologically meaningful representations that generalize across diagnostic and prognostic settings [1, 12, 16, 25].

However, systematic evaluation on realistic, multimodal bedside data remains limited. In continuous monitoring environments, clinicians rely on long-duration single-lead ECG and PPG signals (often spanning minutes to hours in duration), which capture complementary and temporally evolving aspects of cardiovascular physiology [5, 26, 27]. Yet, recent benchmarking studies [2, 3, 18] and FM development papers [16, 21, 23–25, 28] primarily focus on unimodal ECG or PPG data and short signal segments (typically 10 seconds). Moreover, while CSFM [12] demonstrates the promise of large-scale multimodal ECG + PPG pretraining across heterogeneous sensing scenarios, it does not establish a standardized benchmark explicitly designed for controlled comparison of FMs on synchronized long-duration ECG and PPG signals. Consequently, it remains unclear how different types of biosignal FMs perform on realistic multimodal emergency department data, and how these domain-specific FMs compare to general-purpose time-series FMs. In particular, the interaction between multimodal fusion, varying signal lengths, and model scale has not been systematically evaluated, nor has the relative value of learned FM representations versus traditional domain feature baselines been rigorously assessed.

To address these gaps, we introduce *SignalMC-MED*, a clinically grounded benchmark derived from the MC-MED dataset [14, 15], comprising 22,256 adult emergency department visits with 10-minute synchronized ECG and PPG signals (Figure 1). The benchmark includes 20 clinically relevant tasks spanning prediction of demographics, emergency department disposition, laboratory value regression, and detection of *prior* ICD-10 diagnoses, and adopts

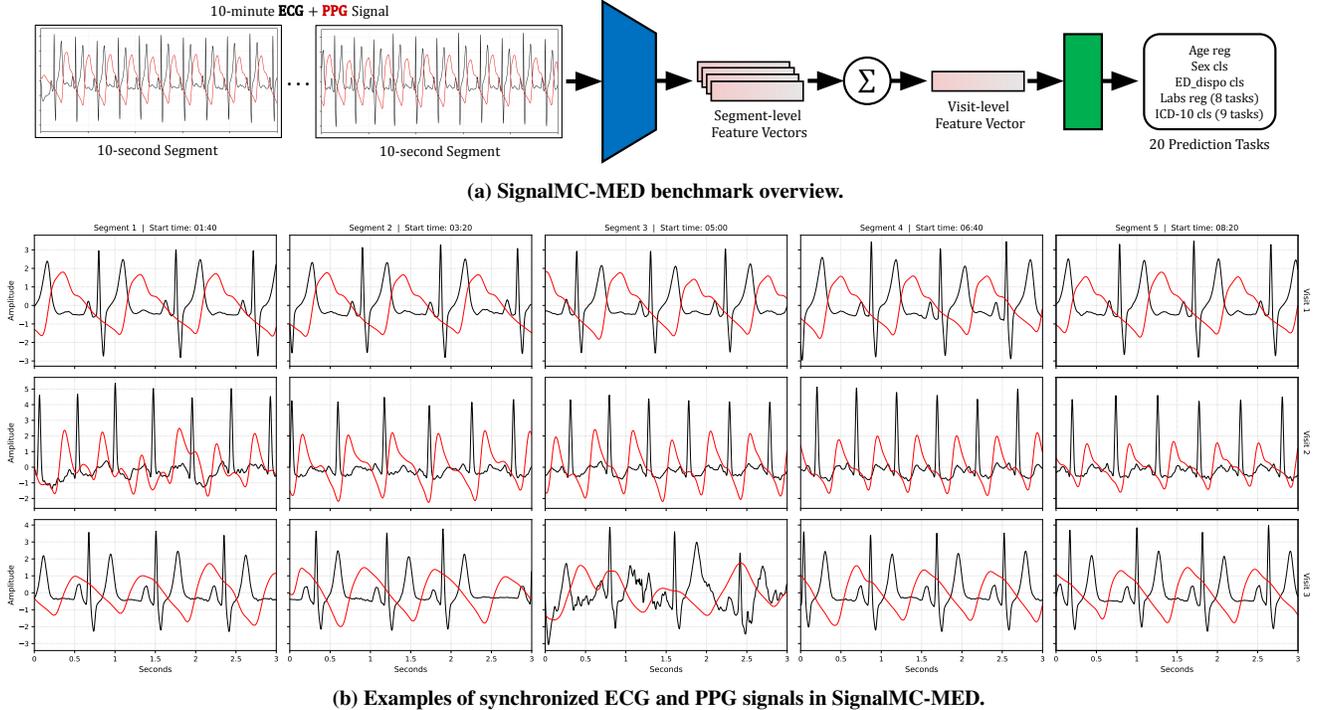


Figure 1. **(a)** Overview of the SignalMC-MED benchmark and evaluation framework. For each of the 22,256 visits, synchronized 10-minute single-lead ECG (**black**) and PPG (**red**) signals are divided into non-overlapping 10-second segments. A frozen FM (**blue**) extracts a feature vector for each 10-second segment. The resulting segment-level features are aggregated via mean pooling (Σ) to form a visit-level representation, which is used to train a linear prediction model (**green**) across 20 downstream tasks: age regression, sex classification, emergency department disposition classification, laboratory value regression (8 tasks), and prior ICD-10 diagnosis classification (9 tasks). **(b)** Examples of synchronized ECG and PPG signals from three different visits. For each visit, five 3-second segments are shown, evenly sampled from the full 10-minute signal (segment start times indicated above each panel). The signals illustrate inter-patient variability in rhythm, morphology, amplitude and noise characteristics. Additional examples are provided in Figure S1 in the supplementary material.

chronological train/val/test splits without patient overlap to approximate real-world deployment conditions. Using this benchmark, we perform a systematic evaluation of representative time-series and biosignal FMs (Table 1).

Our study provides three key contributions: (1) We introduce a large-scale multimodal benchmark for rigorous evaluation of biosignal FMs on synchronized ECG and PPG signals collected in routine emergency care. (2) We conduct a systematic head-to-head evaluation of general time-series FMs, unimodal ECG- and PPG-specific biosignal FMs, a multimodal biosignal FM, and hand-crafted domain features across ECG-only, PPG-only, and ECG + PPG settings. (3) We provide comprehensive empirical evidence clarifying how modality specialization, multimodal fusion, signal duration, model scale, and physiologically grounded domain features influence the performance of biosignal FMs.

Results

We evaluate eight representative approaches spanning two general time-series FMs (MOMENT [10], Chronos-Bolt [4]), four unimodal biosignal FMs (D-BETA [23],

Table 1. **Overview of all evaluated FMs**, including model size (number of parameters) and expected input sampling frequency.

Model Name	Size	Feature Dimension	Sampling Frequency	Pretraining Data
MOMENT-small [10]	35M	512	Mixed	13M time series; 1.2B timestamps
MOMENT-base	110M	768	—	—
MOMENT-large	341M	1024	—	—
Chronos-Bolt-tiny [4]	9M	256	Mixed	100B time series observations
Chronos-Bolt-mini	21M	384	—	—
Chronos-Bolt-small	48M	512	—	—
Chronos-Bolt-base	205M	768	—	—
PaPaGei [24]	6M	512	125Hz	20M 10-sec PPG segments from 13.5K subjects
D-BETA [23]	62M	768	500Hz	0.78M 10-sec, 12-lead ECGs from ~150K subjects
ECGFounder [16]	31M	1024	500Hz	10.8M ECGs (“predominantly 10-sec, 12-lead”) from 1.8M subjects
xECG [18]	57M	1024	100Hz	2.4M 7–10-sec, 12-lead ECGs from ~1.7M subjects + 75 30-min, 12-lead ECGs
CSFM-tiny [12]	43M	768	250Hz	2.3M 10-sec, 12-lead ECGs + 0.27M paired 10-sec PPGs and single-lead ECGs from 1.7M subjects
CSFM-base	109M	768	—	—
CSFM-large	334M	1024	—	—

ECGFounder [16], xECG [18], PaPaGei [24]), one multimodal biosignal FM (CSFM [12]), and hand-crafted ECG and PPG domain features. None of the evaluated FMs were

Table 2. **Main model comparison on the *test* set**, with results aggregated across training using 10%, 25%, 50%, and 100% of the train visits. Regression tasks are evaluated using Pearson correlation (\uparrow), and classification tasks using AUROC (\uparrow). Columns correspond to aggregated performance for: age regression, sex classification, ED disposition classification, laboratory value regression (mean across eight laboratory tests), and prior ICD-10 diagnosis classification (mean across nine diagnosis groups). For each train percentage, training visits are sampled with replacement, and both downstream training and hyperparameter selection are repeated five times. Each table entry represents the mean \pm standard deviation across resampling repetitions, after aggregation across tasks within each category and across train visit percentages. **Bold red** marks the best mean value in each column, **bold** marks the second-best, and underline marks the third-best.

Model Name	Modality	Age reg	Sex cls	ED.dispo cls	Labs reg	ICD-10 cls
Domain Features	ECG-only	0.678 \pm 0.002	0.652 \pm 0.004	0.711 \pm 0.003	0.319 \pm 0.002	0.748 \pm 0.003
Domain Features	PPG-only	0.684 \pm 0.052	0.737 \pm 0.002	0.676 \pm 0.002	0.213 \pm 0.018	0.677 \pm 0.002
Domain Features (concat)	ECG + PPG	0.751 \pm 0.044	0.769 \pm 0.002	0.713 \pm 0.001	0.299 \pm 0.015	0.736 \pm 0.005
MOMENT-base	ECG-only	0.641 \pm 0.002	0.741 \pm 0.002	0.662 \pm 0.004	0.239 \pm 0.005	0.671 \pm 0.005
MOMENT-base	PPG-only	0.702 \pm 0.002	0.724 \pm 0.003	0.676 \pm 0.002	0.240 \pm 0.002	0.670 \pm 0.002
MOMENT-base	ECG + PPG	0.757 \pm 0.002	0.776 \pm 0.001	0.678 \pm 0.003	0.283 \pm 0.004	0.692 \pm 0.004
Chronos-Bolt-small	ECG-only	0.650 \pm 0.004	0.720 \pm 0.002	0.652 \pm 0.004	0.233 \pm 0.002	0.666 \pm 0.004
Chronos-Bolt-small	PPG-only	0.701 \pm 0.003	0.731 \pm 0.002	0.666 \pm 0.002	0.237 \pm 0.002	0.655 \pm 0.004
Chronos-Bolt-small	ECG + PPG	0.745 \pm 0.002	0.765 \pm 0.001	0.669 \pm 0.003	0.270 \pm 0.002	0.682 \pm 0.003
PaPaGei	ECG-only	0.478 \pm 0.003	0.651 \pm 0.002	0.617 \pm 0.002	0.156 \pm 0.001	0.640 \pm 0.003
PaPaGei	PPG-only	0.665 \pm 0.001	0.678 \pm 0.003	0.648 \pm 0.003	0.234 \pm 0.002	0.655 \pm 0.001
PaPaGei	ECG + PPG	0.647 \pm 0.002	0.706 \pm 0.002	0.640 \pm 0.003	0.213 \pm 0.003	0.661 \pm 0.003
D-BETA	ECG-only	0.629 \pm 0.006	0.699 \pm 0.002	0.666 \pm 0.003	0.221 \pm 0.002	0.689 \pm 0.003
D-BETA	PPG-only	0.538 \pm 0.034	0.700 \pm 0.003	0.650 \pm 0.003	0.131 \pm 0.009	0.643 \pm 0.003
D-BETA	ECG + PPG	0.692 \pm 0.010	0.731 \pm 0.002	0.670 \pm 0.004	0.239 \pm 0.005	0.695 \pm 0.005
ECGFounder	ECG-only	0.733 \pm 0.002	0.740 \pm 0.002	0.693 \pm 0.003	0.304 \pm 0.004	0.716 \pm 0.006
ECGFounder	PPG-only	0.681 \pm 0.011	0.724 \pm 0.002	0.673 \pm 0.004	0.229 \pm 0.016	0.667 \pm 0.004
ECGFounder	ECG + PPG	<u>0.799</u> \pm 0.002	<u>0.796</u> \pm 0.001	0.700 \pm 0.004	0.322 \pm 0.002	0.722 \pm 0.005
xECG-10min	ECG-only	0.741 \pm 0.002	0.796 \pm 0.003	0.701 \pm 0.002	0.332 \pm 0.003	0.715 \pm 0.002
xECG-10min	PPG-only	0.753 \pm 0.003	0.751 \pm 0.002	0.692 \pm 0.002	0.258 \pm 0.003	0.681 \pm 0.002
xECG-10min	ECG + PPG	0.803 \pm 0.002	0.822 \pm 0.003	0.707 \pm 0.004	<u>0.348</u> \pm 0.002	0.721 \pm 0.002
CSFM-base	ECG-only	0.752 \pm 0.001	0.791 \pm 0.003	0.709 \pm 0.003	0.350 \pm 0.002	0.729 \pm 0.002
CSFM-base	PPG-only	0.784 \pm 0.002	0.762 \pm 0.002	<u>0.712</u> \pm 0.003	0.303 \pm 0.003	0.701 \pm 0.004
CSFM-base	ECG + PPG	0.818 \pm 0.002	0.823 \pm 0.002	0.721 \pm 0.005	0.374 \pm 0.003	<u>0.735</u> \pm 0.003

pretrained on MC-MED, preventing any overlap between pretraining data and the benchmark. All FMs are evaluated in ECG-only, PPG-only, and ECG + PPG settings, including those pretrained exclusively on a single modality. This unified evaluation allows us to assess cross-modality generalization and determine whether learned representations encode modality-specific morphology or more general cardiovascular dynamics transferable across related biosignals.

All models are used strictly as frozen feature extractors within a unified linear-probing framework to isolate representation quality (see *Methods* and Figure 8). Each 10-minute signal is divided into non-overlapping 10-second segments, segment-level features are extracted, and visit-level representations are obtained via mean aggregation. For shorter-duration experiments (10 seconds to 5 minutes), only the initial contiguous portion of the signal is used. Downstream prediction is performed using ridge regression for regression tasks and regularized logistic regression for classification tasks. For each model, a single regularization hyperparameter per task type is selected on the *val* split.

Performance is reported using Pearson correlation for re-

gression and area under the receiver operating characteristic curve (AUROC) for classification. To assess robustness and data efficiency, downstream models are trained on 10%, 25%, 50%, and 100% of the available *train* visits with five repeated resamplings. Multimodal fusion is implemented via late feature-level fusion by averaging independently extracted ECG-only and PPG-only visit-level representations. The 20 benchmark tasks (Table 4) are grouped into five categories and models are ranked by mean category rank to facilitate cross-model comparison.

Main Model Comparison

Table S1 in the supplementary material reports *val* results for all model variants, with corresponding rankings in Table S2. Based on these, one representative variant each for MOMENT, Chronos-Bolt, xECG, CSFM, and the PPG domain features is selected and evaluated on the *test* set alongside D-BETA, PaPaGei, and ECG domain features (Table 2). Model rankings on *test* are shown separately for ECG-only, PPG-only, and ECG + PPG in Table 3, and jointly across modalities in Table S3. Moreover, Figure 2

shows the *test* performance as a function of the percentage of train visits used, highlighting both absolute performance and data-efficiency trends.

Across modalities, CSFM-base achieves the best overall performance and mean rank, while xECG consistently attains the second-best performance. The two models perform strongly across both regression and classification tasks. In contrast, PaPaGei and D-BETA rank among the weakest performers overall. Notably, PaPaGei, which is the only evaluated PPG-specific FM, ranks last for ECG-only and ECG + PPG, and second-last even for PPG-only. Hand-crafted domain features provide a strong baseline, ranking fourth across all modalities and outperforming multiple FMs. Among ECG-specific models, ECGFounder ranks third for ECG-only and ECG + PPG but drops for PPG-only, whereas xECG, despite ECG-only pretraining, ranks second across modalities and performs strongly even on PPG-only inputs.

Table 4 further reports CSFM-base *val* performance across all 20 benchmark tasks. Performance varies substantially across tasks, with particularly strong results for atrial fibrillation and cardiac device detection, and more modest performance for conditions such as anemia and venous thromboembolism. Results on *val* for all eight evaluated models across the 20 tasks, when training with 10%, 25%, 50%, and 100% of the train visits, are provided in Table S4-S15. Across tasks and train visit percentages, variability across resampling repetitions remains small overall, indicating stable performance estimates and consistent model rankings.

Extended Model Analysis

Beyond the primary benchmark comparison, we further analyze model behavior under different modality settings, signal lengths, and architectural configurations. Figure 3 & S2 show that ECG + PPG fusion consistently improves performance over unimodal inputs for nearly all models and train percentages. Figure 5 & S5 demonstrate that these gains persist across different signal lengths, indicating that multimodal complementarity is not limited to specific temporal resolutions.

Using the full 10-minute signal consistently outperforms shorter segments (Figure 4, S3 & S4), suggesting that longer temporal context provides additional predictive information beyond short-window morphology. In contrast, larger variants within the same model family (e.g., CSFM-tiny/base/large) do *not* reliably outperform their smaller counterparts (Figure S6 & S7), indicating diminishing returns from parameter scaling under frozen linear evaluation.

Figure 6 & S8 evaluate concatenation of hand-crafted domain features with learned FM representations. For ECG, domain-feature augmentation consistently improves perfor-

Table 3. **Main model ranking on the *test* set**, computed separately for (a) ECG-only, (b) PPG-only, and (c) ECG + PPG inputs. Mean model rank (↓) across the five aggregated task categories based on Table 2. A joint ranking across modalities is in Table S3.

Rank	Model Name	Model Ranks	Mean Model Rank (↓)
1	CSFM-base	7,5,5,2,4	4.6
2	xECG-10min	10,3,7,4,8	6.4
3	ECGFounder	11,12,9,7,7	9.2
4	Domain Features	17,23,4,6,1	10.2
5	MOMENT-base	21,11,19,15,16	16.4
6	D-BETA	22,21,18,20,12	18.6
7	Chronos-Bolt-small	19,18,20,18,19	18.8
8	PaPaGei	24,24,24,23,24	23.8

(a) ECG-only.

Rank	Model Name	Model Ranks	Mean Model Rank (↓)
1	CSFM-base	4,9,3,8,9	6.6
2	xECG-10min	6,10,10,12,14	10.4
3	MOMENT-base	12,17,12,13,17	14.2
4	Domain Features	15,13,13,22,15	15.6
5	Chronos-Bolt-small	13,15,17,16,22	16.6
6	ECGFounder	16,16,14,19,18	16.6
7	PaPaGei	18,22,22,17,21	20.0
8	D-BETA	23,20,21,24,23	22.2

(b) PPG-only.

Rank	Model Name	Model Ranks	Mean Model Rank (↓)
1	CSFM-base	1,1,1,1,3	1.4
2	xECG-10min	2,2,6,3,6	3.8
3	ECGFounder	3,4,8,5,5	5.0
4	Domain Features (concat)	8,7,2,9,2	5.6
5	MOMENT-base	5,6,11,10,11	8.6
6	Chronos-Bolt-small	9,8,16,11,13	11.4
7	D-BETA	14,14,15,14,10	13.4
8	PaPaGei	20,19,23,21,20	20.6

(c) ECG + PPG.

mance, whereas improvements are less consistent for PPG. Moreover, simple feature-level fusion slightly outperforms CSFM’s internal multimodal mechanism (Figure S9), and the performance gap between CSFM-base and domain features narrows with longer signals (Figure S10), suggesting that temporal aggregation partially compensates for representational differences.

Finally, UMAP visualizations in Figure 7 & S11 show clearer clustering by age group and atrial fibrillation label for CSFM-base, xECG, and ECGFounder than for general time-series FMs, consistent with their stronger quantitative performance. Although UMAP provides only a qualitative projection, these patterns suggest that domain-specific biosignal FMs can learn more structured and clinically aligned latent representations.

Discussion

Overall, our results demonstrate a clear advantage of domain-specific biosignal FMs over general time-series FMs across ECG-only, PPG-only, and ECG + PPG settings. CSFM-base achieves the best overall performance with xECG consistently second. This suggests that while

Table 4. **SignalMC-MED benchmark task summary and CSFM-base reference performance.** Overview of the 20 included tasks, grouped into two demographics tasks, ED disposition classification, eight laboratory value regression tasks, and nine prior ICD-10 diagnosis classification tasks. CSFM-base performance across modalities on the *val* set when using 100% of the train set visits is shown as a reference for relative task difficulty. Training visits are sampled with replacement, and both downstream training and hyperparameter selection are repeated five times; results are reported as mean \pm standard deviation across resampling repetitions. **Bold red** marks the best mean value in each row. Full results with 100%, 50%, 25%, and 10% of the train visits for all eight evaluated models are provided in Table S4 - S15.

Task	Label / Description	Label Source	ECG-only	PPG-only	ECG + PPG
Age reg	"Age"	visits.csv	0.785 \pm 0.002	0.799 \pm 0.002	0.842 \pm 0.001
Sex cls	"M" vs "F"		0.816 \pm 0.002	0.767 \pm 0.003	0.840 \pm 0.002
ED_dispo cls	"Discharge" vs {"Inpatient", "Observation", "ICU"}	visits.csv	0.744 \pm 0.001	0.726 \pm 0.001	0.750 \pm 0.002
Potassium reg	"POTASSIUM"	labs.csv	0.386 \pm 0.007	0.190 \pm 0.006	0.386 \pm 0.006
Calcium reg	"CALCIUM"		0.209 \pm 0.010	0.183 \pm 0.015	0.265 \pm 0.010
eGFR reg	"EGFR REFIT WITHOUT RACE (2021)"		0.591 \pm 0.006	0.567 \pm 0.005	0.613 \pm 0.003
Glucose reg	"GLUCOSE"		0.262 \pm 0.017	0.297 \pm 0.009	0.322 \pm 0.010
Hemoglobin reg	"HEMOGLOBIN (HGB)"		0.373 \pm 0.007	0.392 \pm 0.009	0.423 \pm 0.007
Albumin reg	"ALBUMIN"		0.428 \pm 0.004	0.394 \pm 0.006	0.450 \pm 0.006
BUN reg	"BLOOD UREA NITROGEN (BUN)"		0.394 \pm 0.005	0.368 \pm 0.004	0.410 \pm 0.004
Sodium reg	"SODIUM"		0.260 \pm 0.009	0.227 \pm 0.012	0.272 \pm 0.006
I48 cls	Atrial Fibrillation Detection	pmh.csv	0.861 \pm 0.005	0.814 \pm 0.004	0.857 \pm 0.005
I50 cls	Heart Failure Detection		0.840 \pm 0.004	0.773 \pm 0.006	0.839 \pm 0.004
{Z95.0, Z95.2, Z95.3} cls	Cardiac Device Detection		0.854 \pm 0.010	0.735 \pm 0.016	0.852 \pm 0.009
N18 cls	Chronic Kidney Disease Detection		0.783 \pm 0.005	0.748 \pm 0.004	0.782 \pm 0.003
{E10, E11} cls	Diabetes Mellitus Detection		0.696 \pm 0.005	0.702 \pm 0.012	0.720 \pm 0.006
{G47.3, E66.2} cls	Sleep-Related Breathing Disorder Detection		0.686 \pm 0.016	0.715 \pm 0.006	0.715 \pm 0.017
{D50, D51, D52, D53, D64} cls	Anemia Detection		0.671 \pm 0.010	0.665 \pm 0.009	0.669 \pm 0.011
{J96, R09.02, R06.0} cls	Respiratory Compromise Detection		0.676 \pm 0.004	0.673 \pm 0.004	0.685 \pm 0.002
{I26, I80, I82} cls	Venous Thromboembolism Detection		0.675 \pm 0.011	0.659 \pm 0.011	0.685 \pm 0.015

general time-series FMs may learn broadly transferable sequence patterns, they lack inductive biases tailored to cardiovascular physiology. However, the poor performance of PaPaGei even in the PPG-only setting shows that modality-specific pretraining alone does not guarantee strong performance. Notably, although PaPaGei is pretrained on 20 million PPG segments, these originate from only 13,500 subjects, whereas CSFM, ECGFounder, and xECG are trained on data from more than 1.5 million individuals (Table 1). The substantially greater subject-level diversity in the latter models may expose them to a wider range of physiological variability and signal conditions, which likely contributes to improved downstream generalization beyond what can be achieved through segment-level scale alone.

Multimodal ECG + PPG fusion yields robust and consistent gains across tasks, train percentages, and signal lengths, including for models pretrained exclusively on ECG, highlighting that PPG provides complementary physiological information beyond what can be captured from ECG alone. Notably, the simple late-fusion strategy of averaging unimodal features slightly outperforms CSFM’s internal multimodal mechanism. This may reflect several factors, including the relative scarcity of paired ECG + PPG segments compared to ECG-only data in CSFM’s multimodal pretraining (Table 1), and the absence of explicit cross-modality alignment losses. Under frozen linear evaluation, preserving unimodal structure and combining it at the feature level may also allow the downstream model to exploit complementary information more flexibly than relying on

a fixed internal fusion mechanism. These findings suggest that effective multimodal integration remains a non-trivial design challenge.

Signal duration proves more impactful than architectural scale under our frozen linear evaluation: extending inputs to the full 10-minute recordings produces larger and more reliable improvements than increasing parameter count, and larger model variants do not consistently outperform mid-sized models. When representations are projected to a 256-dimensional space using PCA, performance trends remain similar (Figure S6 & S7), indicating that the occasional underperformance of larger model variants is not attributable to differences in feature dimensionality alone. This scaling saturation may instead reflect a representational bottleneck under linear probing: once key rhythm and morphology features are linearly accessible, additional capacity may encode information that is not exploitable without fine-tuning. It may also indicate that model size is not the primary constraint at current pretraining data scales.

Hand-crafted domain features provide a strong baseline and usually outperform both MOMENT and Chronos-Bolt, highlighting the continued value of physiologically grounded feature engineering. Importantly, feature concatenation experiments further reveal complementary strengths between learned and hand-crafted representations. For ECG, augmenting FM embeddings with domain features consistently improves performance, even for models that already outperform standalone domain features, suggesting that explicit morphological descriptors add informa-

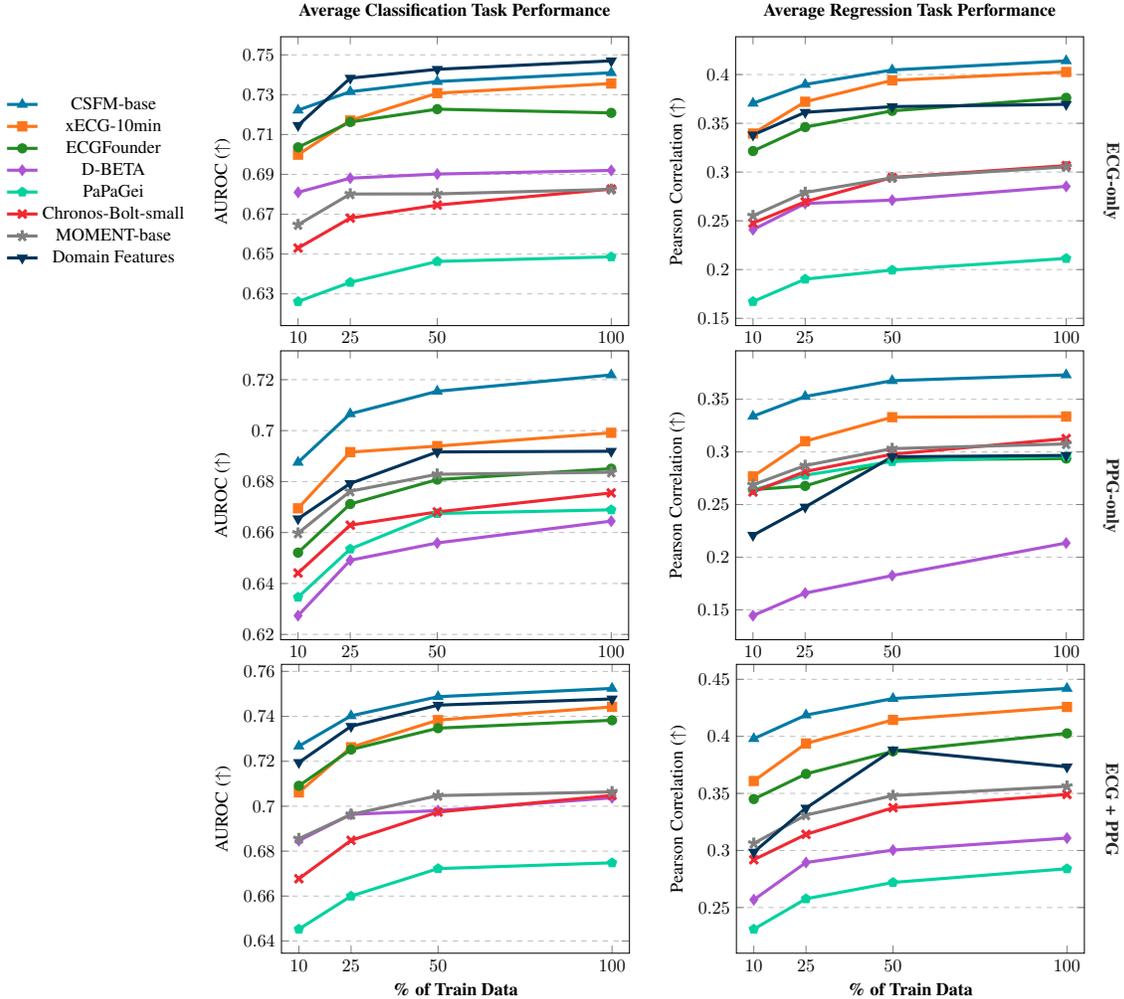


Figure 2. **Main model comparison on the test set**, with average classification task performance (eleven tasks, left column) and regression task performance (nine tasks, right column) shown as a function of the percentage of train visits used (10%, 25%, 50%, 100%). Rows correspond to ECG-only (top), PPG-only (middle), and ECG + PPG (bottom) inputs. Each curve represents one evaluated model and shows the mean performance across five resampling repetitions (standard deviations omitted for visual clarity).

tion beyond learned FM representations. The comparatively weaker and less consistent gains for PPG indicate that our ECG-derived features encode particularly informative structure. These findings suggest that current ECG FMs may not fully internalize certain established morphological abstractions, and that incorporating auxiliary objectives or physiologically informed inductive biases during pretraining could help ensure that learned FM representations subsume strong handcrafted features rather than merely complement them.

Interestingly, xECG performs strongly even on PPG-only inputs despite ECG-only pretraining, suggesting partial transferability of learned temporal and rhythm-level representations across related cardiovascular signals. Both ECG and PPG reflect manifestations of the same underlying cardiac cycle, and models trained to capture long-

range temporal structure and inter-beat dynamics in ECG may therefore learn higher-level cardiovascular representations that remain informative for PPG. In contrast, the PPG-specific model PaPaGei does not generalize well beyond its training modality, which may in part reflect the limited subject-level diversity in its pretraining data. Together, these findings suggest that large-scale exposure to diverse cardiovascular dynamics may be more important for cross-modality robustness than strict modality specialization alone.

Finally, variability in classification task performance does not follow class prevalence alone (Table S17). For example, cardiac device detection (2.4% test prevalence) and atrial fibrillation detection (7.6%) achieve among the highest AUROC values, whereas more prevalent conditions such as anemia (11.1%) and respiratory compromise (13.1%)

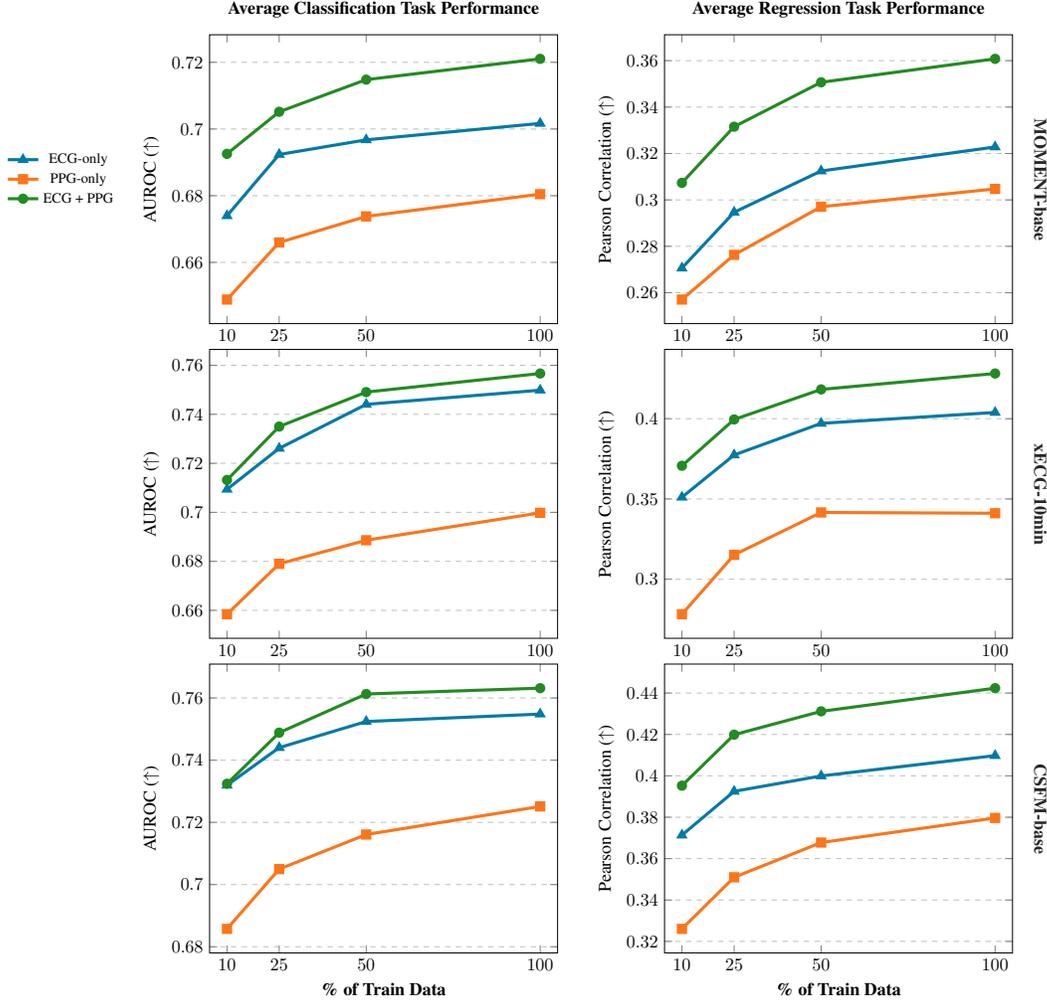


Figure 3. **Multimodal fusion comparison on the val set**, comparing the ECG-only, PPG-only, and ECG + PPG performance for three representative models (MOMENT-base, xECG-10min, CSFM-base). Average classification (left) and regression (right) task performance is shown as a function of the percentage of train visits used. Results for the five remaining models are provided in Figure S2.

show more modest performance. This indicates that predictive accuracy is driven primarily by the strength and specificity of physiological signatures present in ECG and PPG signals rather than by label frequency. Indeed, conditions such as atrial fibrillation and cardiac device implantation are known to produce clear and direct waveform-level manifestations, whereas anemia does not have a specific electrocardiographic signature. Venous thromboembolism may also occur without overt cardiac strain, consistent with the more modest performance observed for this task.

This study has several limitations. All evaluations rely on frozen representations with linear downstream models, and results may differ under end-to-end fine-tuning or non-linear adaptation. The benchmark focuses on the first 10 minutes of ECG and PPG recordings and employs simple temporal aggregation strategies. Future work should ex-

tend the evaluation to continuous monitoring scenarios with hours or days of data per patient to better assess long-term temporal modeling. We evaluate seven representative FMs spanning general time-series models and both unimodal and multimodal biosignal models, but many additional models could be included in future comparisons. Although chronological splits prevent temporal leakage and patient overlap, all data originate from a single healthcare system. Generalizability across institutions, acquisition devices, and patient populations therefore remains to be evaluated. Finally, while the task suite spans 20 clinically relevant outcomes, it does not cover the full spectrum of cardiovascular and systemic conditions, nor does it assess prospective clinical integration. Future work should also explore more advanced multimodal fusion strategies and temporal aggregation methods, and evaluate robustness under distribution

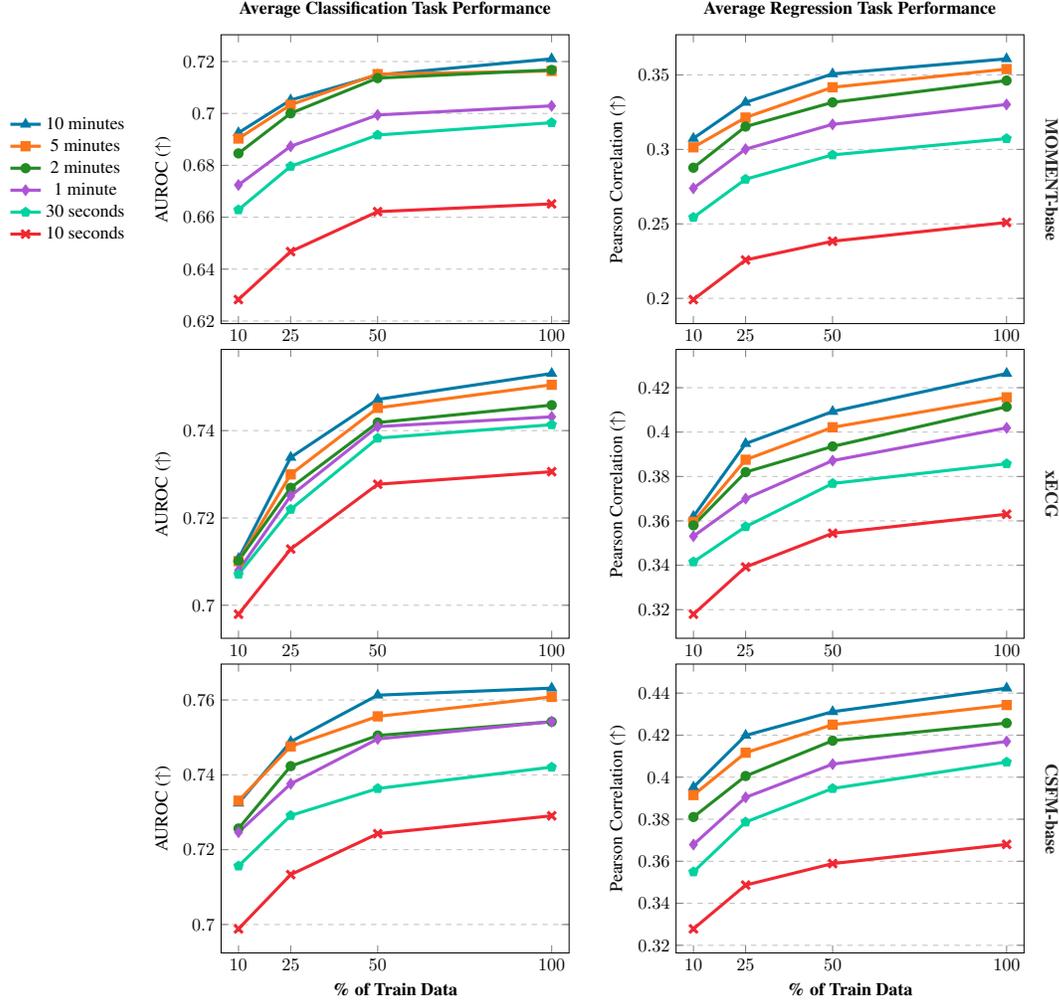


Figure 4. **Signal length comparison on the val set**, showing how the *ECG + PPG* performance is affected by increasing the signal length from 10 seconds to 10 minutes for three representative models (MOMENT-base, xECG, CSFM-base). Results for the five remaining models are provided in Figure S3, and for ECG and PPG domain features in Figure S4.

shifts and real-world deployment settings.

The main actionable takeaways from our study are: (1) Domain-specific biosignal FMs should be preferred over general time-series FMs for cardiovascular prediction tasks. (2) Multimodal ECG and PPG data should be leveraged whenever available, as even simple late-fusion strategies yield consistent improvements over unimodal inputs. (3) Longer signal segments should be prioritized, as increased temporal coverage (up to the full 10-minute signal) consistently improves performance across models and tasks. (4) Larger FM variants should not be assumed to outperform smaller ones without task-specific validation. (5) Hand-crafted ECG domain features should be included as a simple yet strong baseline to verify that FMs provide meaningful gains, and can offer complementary benefits when combined with learned FM representations.

Methods

We describe the construction of the SignalMC-MED benchmark, the downstream prediction tasks, and the unified evaluation framework used to compare FMs. As illustrated in Figure 8, each visit comprises a synchronized 10-minute single-lead ECG and PPG signal that is segmented, encoded using frozen FMs, aggregated into visit-level representations, and evaluated using linear models across 20 clinically relevant tasks. The following sections detail the dataset construction, preprocessing, task definitions, evaluation protocol, and the evaluated FMs and baselines.

Dataset Construction

SignalMC-MED is derived from the MC-MED dataset [14, 15], which consists of data from 118,385 emergency de-

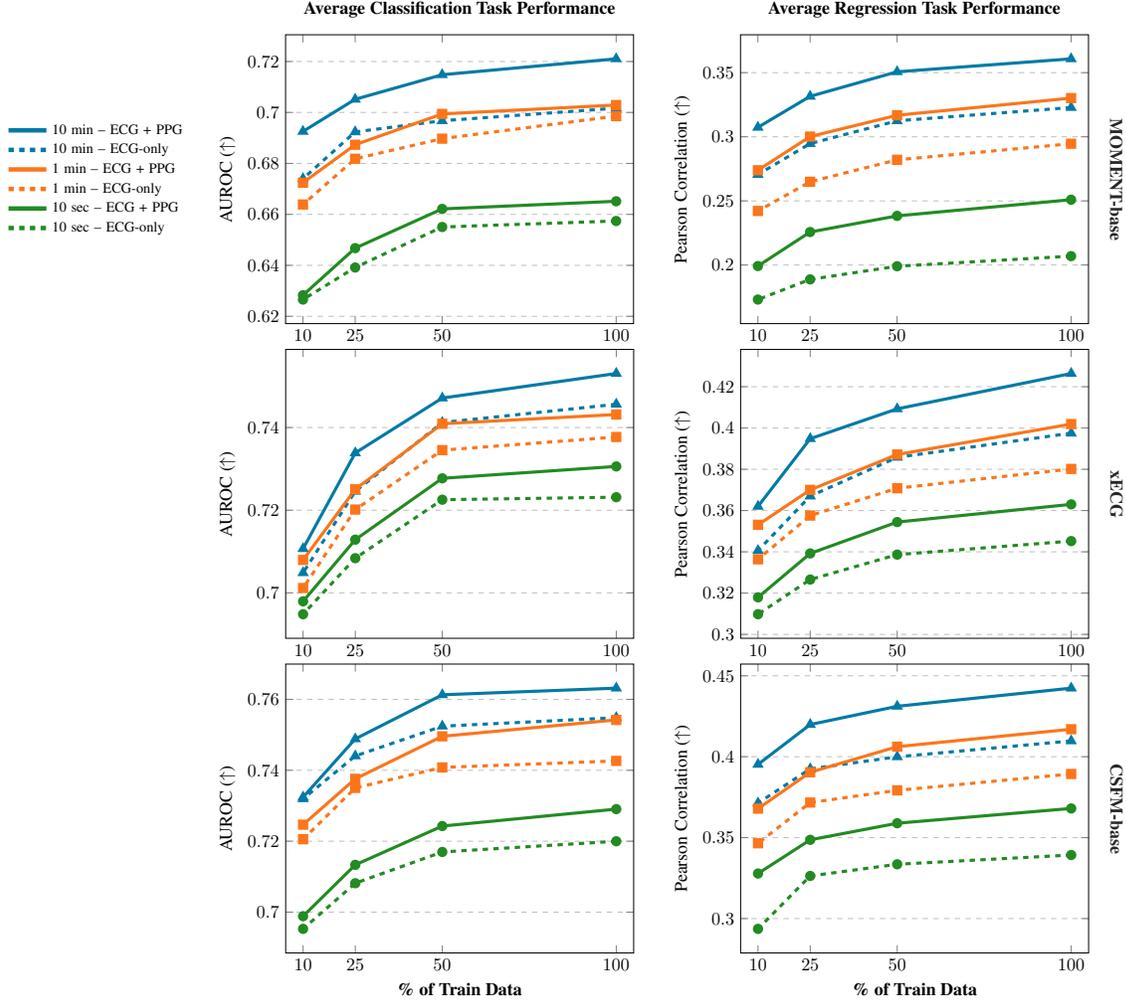


Figure 5. **Multimodal fusion comparison across signal lengths on the *val* set**, comparing the ECG-only (dashed lines) and ECG + PPG (solid lines) performance with multiple signal lengths (10 seconds, 1 minute, 10 minutes) for three representative models (MOMENT-base, xECG, CSFM-base). Results for the five remaining models are provided in Figure S5.

partment visits involving 70,545 unique adult patients, collected at Stanford Health Care between 2020 and 2022. MC-MED includes continuously acquired single-lead ECG and PPG recordings captured using bedside patient monitors, together with structured clinical information such as visit outcomes, vital signs, prior medical history, and laboratory test results. Among the 54,024 visits with recordings of both ECG and PPG, we identify 24,567 visits with at least 10 minutes of temporally overlapping ECG and PPG data within the first recorded signal segment. After signal preprocessing and quality control, this yields 22,256 visits with valid synchronized ECG and PPG signals, corresponding to 17,841 unique patients. For each visit, a contiguous 10-minute ECG/PPG segment is extracted and used for all downstream analyses.

We adopt the chronological train/val/test splits provided

by MC-MED to reflect a realistic clinical deployment scenario. All *val* visits occur after the final *train* visit, and all *test* visits occur after the final *val* visit, with no patient appearing in more than one split. This gives 15,832 *train* visits (13,011 unique patients), 2,030 *val* visits (1,882 unique patients), and 1,995 *test* visits (1,870 unique patients).

Signal Preprocessing & Segmentation Each 10-minute ECG and PPG signal undergoes standardized preprocessing consisting of resampling to a common frequency, signal cleaning using NeuroKit2 [19], and Z-score normalization. ECG and PPG signals are resampled to the same frequency, matched to the pretraining sampling frequency expected by each model (Table 1), with 250 Hz used as the default where applicable. The preprocessed 10-minute signals are then divided into non-overlapping 10-second segments, which form the basic input units for all FMs. For experiments in-

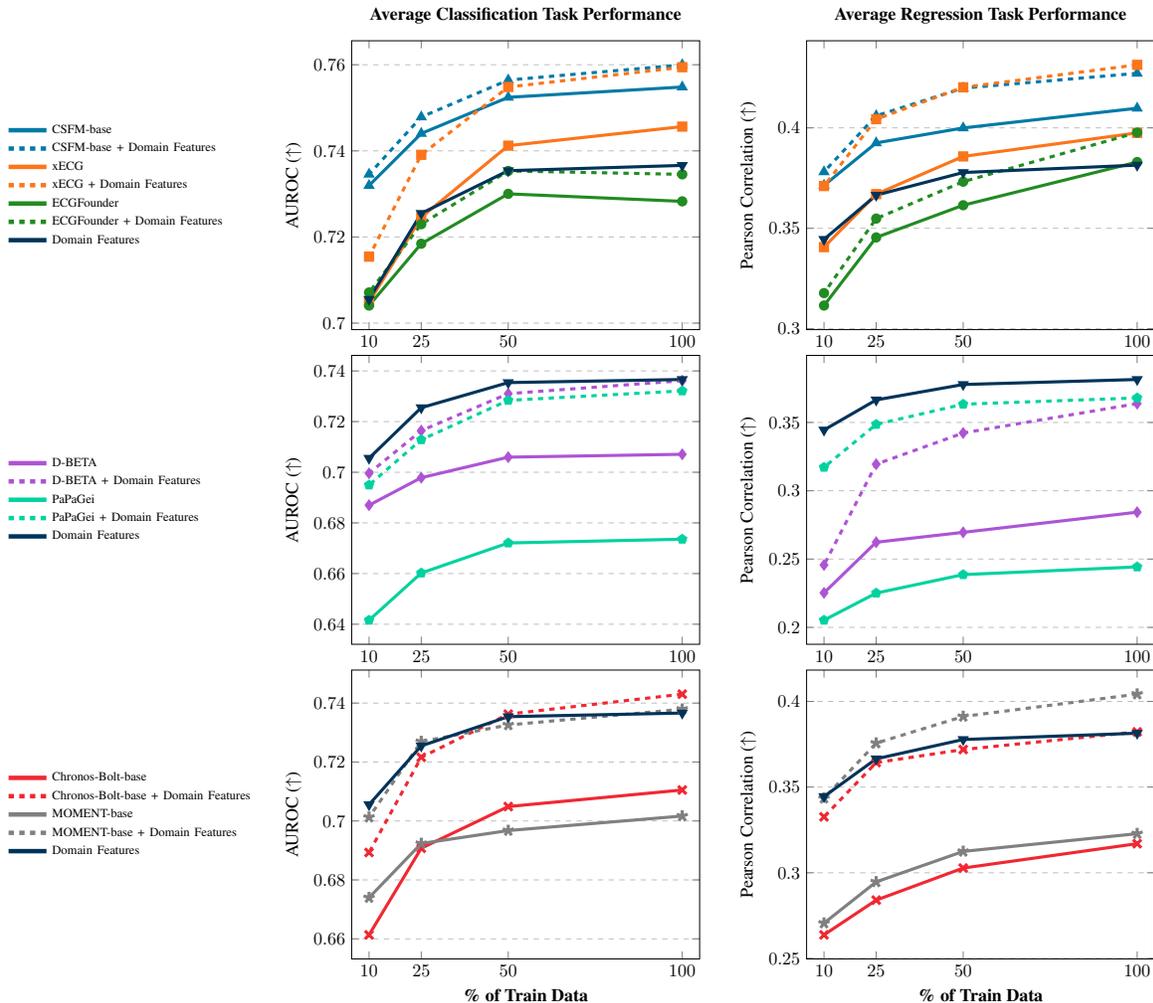


Figure 6. **Effect of concatenating ECG domain features with learned FM representations on the *val* set.** Comparison of ECG-only performance before (solid lines) and after (dashed lines) concatenation of hand-crafted ECG domain features to the frozen FM feature vectors, shown for each of the seven evaluated FMs. The same comparison for PPG domain features is provided in Figure S8.

volving shorter signal durations, only the initial contiguous subset of segments corresponding to the desired duration is used. As signal quality control, hand-crafted domain features are also extracted from both the ECG and PPG signals (see *Domain Features* below for details). Visits for which any preprocessing step fails are excluded from the analysis.

Benchmark Tasks

SignalMC-MED includes 20 downstream prediction tasks spanning both regression and binary classification. Regression tasks include patient age and eight laboratory measurements obtained during the emergency department visit. Classification tasks include patient sex, emergency department disposition, and the presence of specific ICD-10 diagnoses in the patient’s prior medical history, grouped into nine clinically meaningful code sets. These tasks were se-

lected to capture a diverse range of physiological, demographic, and clinically relevant outcomes. Table 4 summarizes all 20 benchmark tasks and reports CSFM-base performance across them. The positive class prevalence for each of the binary classification tasks across the *train*, *val*, and *test* splits is reported in Table S17. These results provide a reference for relative task difficulty and highlight the heterogeneity of signal predictability across clinical outcomes.

Demographics We include two demographic tasks: age regression and sex classification. Age is treated as a continuous target representing patient age in years at the time of the visit, and sex is formulated as a binary classification task using the recorded biological sex. Labels for both tasks are extracted from the `visits.csv` file in MC-MED.

ED Disposition Emergency department disposition is formulated as a binary classification task distinguishing vis-

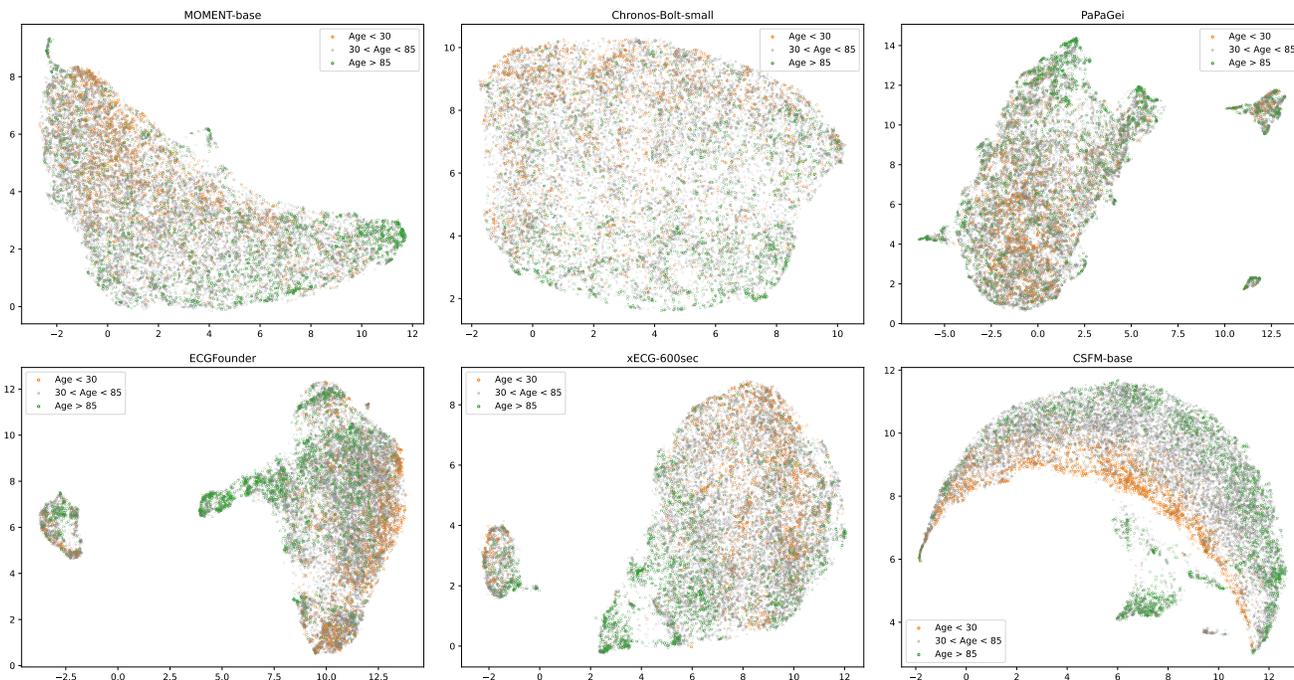


Figure 7. **UMAP visualizations colored by age group.** UMAP [20] projections of the *train* set *ECG + PPG* feature vectors for six selected models, colored according to three age groups (< 30 years, 30 – 85 years, > 85 years). CSFM-base, xECG-10min and ECGFounder exhibit clearer age-related structure compared to general time-series FMs (MOMENT-base, Chronos-Bolt-small) and PaPaGei.

its resulting in discharge from those resulting in inpatient admission, observation status, or intensive care unit admission. Disposition labels are extracted from `visits.csv`.

Laboratory Measurements Eight laboratory regression tasks corresponding to commonly ordered laboratory tests are included (e.g., electrolytes, renal function markers, and metabolic indicators). For each task, the laboratory value measured during the visit is used as the regression target, extracted from `labs.csv`. To ensure unambiguous supervision, only visits with exactly one recorded measurement of the corresponding laboratory test are included. The resulting number of *train*, *val*, and *test* visits therefore varies slightly by task, as reported in Table S16.

ICD-10 Codes Nine binary classification tasks are defined for detecting whether a patient has a documented history of specific ICD-10 diagnoses *prior to* the emergency department visit. Diagnoses are grouped into clinically meaningful categories, including atrial fibrillation, heart failure, cardiac device implantation, chronic kidney disease, diabetes mellitus, sleep-related breathing disorder, anemia, respiratory compromise, and venous thromboembolism. Labels are derived from the patient prior medical history in `pmh.csv`.

Evaluation Framework

All FMs are evaluated using a unified linear probing framework designed to isolate representation quality, as illustrated in Figure 8. Models are used exclusively as frozen feature extractors, and no fine-tuning is performed on downstream tasks.

Feature Extraction & Aggregation Each 10-minute signal is divided into non-overlapping 10-second segments. For each segment, a fixed-dimensional feature vector is extracted from the FM. Segment-level features are aggregated by computing the mean feature vector across all segments within the selected signal duration (i.e., averaging all 60 segment feature vectors for the full 10-minute signal, or averaging the first 6 feature vectors for a 1-minute signal). For the multimodal setting, simple feature-level fusion is performed by averaging the ECG-only and PPG-only visit-level feature vectors. For each visit and FM, this yields an ECG-only, a PPG-only, and an ECG + PPG feature vector of the same dimensionality.

Downstream Prediction Models Downstream prediction is performed using linear models on top of the aggregated visit-level feature vectors. Ridge regression is used for regression tasks and logistic regression for classification tasks. To minimize task-specific tuning, a single regularization hyperparameter is selected per FM. The ridge penalty is

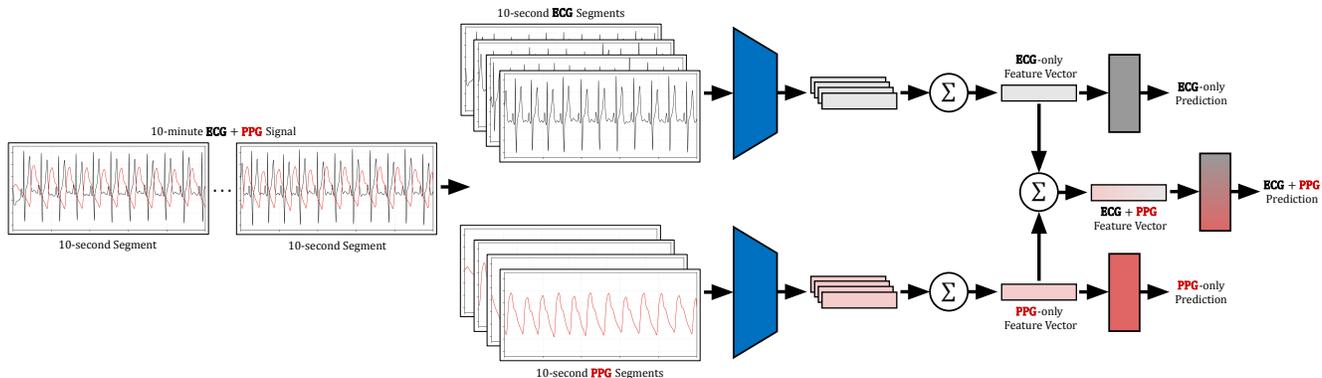


Figure 8. **Detailed evaluation framework for SignalMC-MED.** A 10-minute synchronized single-lead ECG (black) and PPG (red) signal is divided into non-overlapping 10-second segments. ECG and PPG segments are processed separately by a frozen FM (blue) to extract a feature vector for each 10-second segment. The evaluated FM is applied to both ECG and PPG inputs, regardless of its pretraining modality. The resulting segment-level features are aggregated via mean pooling (Σ) to form ECG-only and PPG-only visit-level feature vectors, which are then averaged to obtain an ECG + PPG feature vector of the same dimensionality. Separate linear models are trained on the ECG-only, PPG-only, and ECG + PPG representations to generate task-specific predictions across the 20 benchmark tasks.

optimized on the *val* split using the age regression task, and the logistic regression inverse regularization strength is optimized using the sex classification task. These hyperparameters are then fixed and applied across all remaining tasks. Hyperparameter selection is performed separately for ECG-only, PPG-only, and ECG + PPG representations, as well as for each evaluated signal duration. Apart from this regularization selection, all downstream training settings are held constant across modalities, signal lengths, and FMs to ensure fair comparison. The implementation uses scikit-learn [22] (version 1.7.2) with default settings. Feature vectors are standardized using the mean and standard deviation computed on the *train* split, and the same transformation is applied to both *val* and *test*.

Evaluation Metrics & Reporting Model performance is evaluated using Pearson correlation for regression tasks and AUROC for classification tasks. To assess data efficiency and robustness, downstream models are trained using 10%, 25%, 50%, and 100% of the available *train* visits. For each percentage, training visits are sampled with replacement from the *train* split (including the 100% setting), and both downstream training and hyperparameter selection are repeated five times. Results across all 20 benchmark tasks (Table 4, Table S4 - S15) are reported as mean \pm standard deviation over the five resampling repetitions, for a given train visit percentage.

For summary comparisons (Table 2), the 20 benchmark tasks are grouped into five categories: age regression, sex classification, emergency department disposition classification, laboratory value regression (aggregated across eight laboratory tests), and ICD-10 diagnosis classification (aggregated across nine diagnosis groups). For each resampling repetition, performance is first aggregated across tasks

within each category and across train visit percentages, yielding one aggregate value per category. The reported results correspond to the mean \pm standard deviation across the five resampling repetitions. This produces five aggregate mean performance values per FM and modality. Based on these, models are ranked separately for ECG-only, PPG-only, and ECG + PPG by computing ranks within each of the five task categories and then averaging these ranks across categories (Table 3). The mean model ranks are used to facilitate cross-model comparison while reducing sensitivity to any individual task.

For figure-based comparisons (Figure 2 - 6, Figure S2 - S10), classification performance is aggregated across the eleven binary classification tasks and regression performance across the nine regression tasks. For clarity of visualization, figures display only the mean performance across the five resampling repetitions, omitting standard deviation bars to avoid visual clutter.

Evaluated Models

We evaluate eight representative approaches spanning general time-series FMs, ECG-specific and PPG-specific biosignal FMs, a multimodal ECG + PPG FM, and hand-crafted domain features. An overview of the evaluated FMs is provided in Table 1, including model size, pretraining data, feature dimensionality, and expected input sampling frequency. None of the evaluated FMs were pretrained on MC-MED, ensuring strict separation between model pretraining and benchmark evaluation. When multiple model sizes are available within a family, the variant used for primary comparisons is selected based on validation performance (Table S1).

General Time-Series FMs We evaluate MOMENT [10] and Chronos-Bolt [4], which are modality-agnostic time-series FMs pretrained on large-scale heterogeneous time-series corpora. Multiple model sizes are considered (Table 1). Because these models were pretrained on data with varying sampling frequencies and do not specify an expected input frequency, we resample ECG and PPG signals to the default 250 Hz across all experiments. For Chronos-Bolt, segment-level feature vectors are extracted by using the *ChronosBoltPipeline.embed* function and mean pooling across all patch token embeddings. For MOMENT, we use the *MOMENTPipeline* in its default embedding configuration.

Unimodal Biosignal FMs We evaluate four unimodal biosignal FMs: three ECG-specific models (D-BETA [23], ECGFounder [16] and xECG [18]) and one PPG-specific model (PaPaGei [24]). The ECG-specific models are pretrained on large-scale ECG datasets and incorporate architectural and training strategies tailored to electrocardiographic morphology and rhythm modeling. PaPaGei is pretrained on large-scale PPG datasets comprising millions of short PPG segments. We use the PaPaGei-S variant across all experiments. Inputs are resampled to match each model’s expected sampling frequency according to Table 1. Because xECG is based on the xLSTM architecture [6], which supports long input sequences, we additionally evaluate a variant that directly processes the entire 10-minute signal and extracts a single visit-level feature vector. This long-input variant is denoted *xECG-10min*. Although all these models are pretrained exclusively on single-modality data, we evaluate them on both ECG-only, PPG-only, and ECG + PPG inputs to assess cross-modality generalization.

For D-BETA, segment-level feature vectors are extracted using the *get_ecg_feats* function, where the ECG/PPG segment is provided as the second channel and the remaining 11 channels are set to zero. For ECGFounder, we use the single-lead model version and extract features using *ft_1lead_ECGFounder* with *return_features* set to true. For both xECG and xECG-10min, the ECG/PPG segment is provided as the second channel with the remaining leads set to zero. For PaPaGei, we use *ResNet1DMoE* configured with a single input channel.

Multimodal Biosignal FM We evaluate CSFM [12], a multimodal biosignal FM pretrained on more than 2 million ECGs and over 250,000 paired ECG and PPG recordings. CSFM is designed to learn both shared and modality-specific representations across cardiovascular signals. Inputs are resampled to the 250 Hz sampling frequency used during pretraining, and multiple model sizes are considered (Table 1). CSFM supports ECG-only, PPG-only, and joint ECG + PPG inputs via an internal multimodal fusion mechanism. For consistency with the other evaluated

models, primary results are reported using simple feature-level fusion obtained by averaging the unimodal representations. A direct comparison between this feature-level fusion and CSFM’s internal fusion mechanism is provided in Figure S9. For ECG-only and PPG-only inputs, the respective segment is provided as a single-channel input with $\text{channel} = [1]$ (ECG) or $\text{channel} = [12]$ (PPG). For joint ECG + PPG inputs, the two segments are concatenated and provided as a two-channel input with $\text{channel} = [1, 12]$.

Domain Features As a simple yet strong baseline, we evaluate hand-crafted domain features extracted independently from ECG and PPG signals. Following CSFM [12], we compute a 54-dimensional ECG feature vector using NeuroKit2 [19] and a 306-dimensional PPG feature vector using pyPPG [8]. Signals are resampled to 250 Hz prior to feature extraction for consistency. Features are extracted at the segment level and aggregated using the same mean pooling strategy applied to FM representations. Standard 10-second segments are used for ECG, whereas 20-second segments are used for PPG to satisfy pyPPG’s minimum input length requirement. For PPG, we also consider 60-second and 120-second segments (Table S1). In the multimodal ECG + PPG setting, ECG and PPG domain feature vectors are concatenated to form a combined 360-dimensional representation.

Ethics Statement

This study exclusively utilizes publicly available, de-identified data from the MC-MED dataset [14, 15], and complies with the data use policies governing MC-MED.

Data Availability

SignalMC-MED is derived from the MC-MED dataset which is available on PhysioNet [9] at <https://physionet.org/content/mc-med>.

Code Availability

The code for preprocessing, feature extraction and model evaluation is available at <https://github.com/fregu856/SignalMC-MED>. Further implementation details are available from FKG upon reasonable request.

Acknowledgments

DAC was funded by an NIHR Research Professorship; a Royal Academy of Engineering Research Chair; and the InnoHK Hong Kong Centre for Cerebro-cardiovascular Engineering (COCHE); and was supported by the National Institute for Health and Care Research (NIHR) Oxford Biomedical Research Centre (BRC) and the Pandemic Sciences Institute at the University of Oxford. XG was supported by funding from the Engineering and Physical Sciences Research Council IAA (EP/X525777/1), the Medical Research

Council (MR/X50273X/1), and a Smart Data Research UK Fellowship (UKRI4004). This study was also supported by the NIHR Health Protection Research Unit in Healthcare Associated Infections and Antimicrobial Resistance (NIHR207397), a partnership between the UK Health Security Agency (UKHSA), the University of Oxford, and the NIHR Oxford Biomedical Research Centre. The views expressed are those of the authors and not necessarily those of the NIHR, UKHSA or the Department of Health and Social Care.

Author Contributions

FKG was responsible for project conceptualization, software implementation, preparation of figures and tables, and manuscript drafting. XG, MC, and PP contributed to the design of experiments and to the interpretation and presentation of results. DAC and XG supervised the project. DAC and DWE acquired funding. All authors contributed to manuscript revision and finalization.

Competing Interests

The authors declare no competing interests.

Declaration of Generative AI Use

During the preparation of this manuscript, the authors used ChatGPT 5.2 to assist with language editing and drafting. All content produced using this tool was critically reviewed, edited and validated by the authors, who take full responsibility for the final content of the manuscript.

References

- [1] Salar Abbaspourzad, Oussama Elachqar, Andrew Miller, Saba Emrani, Udhayakumar Nallasamy, and Ian Shapiro. Large-scale training of foundation models for wearable biosignals. In *International Conference on Learning Representations (ICLR)*, 2024. [1](#)
- [2] M A Al-Masud, Juan Lopez Alcaraz, and Nils Strodthoff. Benchmarking ECG foundational models: A reality check across clinical tasks. In *International Conference on Learning Representations (ICLR)*, 2026. [1](#)
- [3] Juan Miguel Lopez Alcaraz, Hjalmar Bouma, and Nils Strodthoff. Enhancing clinical decision support with physiological waveforms—a multimodal benchmark in emergency care. *Computers in Biology and Medicine*, 192:110196, 2025. [1](#)
- [4] Abdul Fatir Ansari, Lorenzo Stella, Ali Caner Turkmen, Xiyuan Zhang, Pedro Mercado, Huibin Shen, Oleksandr Shchur, Syama Sundar Rangapuram, Sebastian Pineda Arango, Shubham Kapoor, Jasper Zschiegner, Danielle C. Maddix, Hao Wang, Michael W. Mahoney, Kari Torkkola, Andrew Gordon Wilson, Michael Bohlke-Schneider, and Bernie Wang. Chronos: Learning the language of time series. *Transactions on Machine Learning Research (TMLR)*, 2024. [2](#), [13](#)
- [5] Karim Bayoumy, Mohammed Gaber, Abdallah Elshafeey, Omar Mhaimeed, Elizabeth H Dineen, Francoise A Marvel, Seth S Martin, Evan D Muse, Mintu P Turakhia, Khaldoun G Tarakji, and Mohamed B Elshazly. Smart wearable devices in cardiovascular care: where we are and how to move forward. *Nature Reviews Cardiology*, 18(8):581–599, 2021. [1](#)
- [6] Maximilian Beck, Korbinian Pöppel, Markus Spanning, Andreas Auer, Oleksandra Prudnikova, Michael Kopp, Günter Klambauer, Johannes Brandstetter, and Sepp Hochreiter. xLSTM: Extended long short-term memory. *Advances in Neural Information Processing Systems (NeurIPS)*, 37:107547–107603, 2024. [13](#)
- [7] Rishi Bommasani, Drew A Hudson, Ehsan Adeli, Russ Altman, Simran Arora, Sydney von Arx, Michael S Bernstein, Jeannette Bohg, Antoine Bosselut, Emma Brunskill, et al. On the opportunities and risks of foundation models. *arXiv preprint arXiv:2108.07258*, 2021. [1](#)
- [8] Márton Á Goda, Peter H Charlton, and Joachim A Behar. pyPPG: a python toolbox for comprehensive photoplethysmography signal analysis. *Physiological Measurement*, 45(4):045001, 2024. [13](#)
- [9] Ary L Goldberger, Luis AN Amaral, Leon Glass, Jeffrey M Hausdorff, Plamen Ch Ivanov, Roger G Mark, Joseph E Mietus, George B Moody, Chung-Kang Peng, and H Eugene Stanley. PhysioBank, PhysioToolkit, and PhysioNet: Components of a new research resource for complex physiological signals. *Circulation*, 101(23):e215–e220, 2000. [13](#)
- [10] Mononito Goswami, Konrad Szafer, Arjun Choudhry, Yifu Cai, Shuo Li, and Artur Dubrawski. MOMENT: A family of open time-series foundation models. In *International Conference on Machine Learning (ICML)*, pages 16115–16152, 2024. [2](#), [13](#)
- [11] Xiao Gu, Yuxuan Shu, Jinpei Han, Yuxuan Liu, Zhang-daihong Liu, James Anibal, Veer Sangha, Edward Phillips, Bradley Segal, Hang Yuan, et al. Foundation models for biosignals: A survey. *Authorea Preprints*, 2025. [1](#)
- [12] Xiao Gu, Wei Tang, Jinpei Han, Veer Sangha, Fenglin Liu, Shreyank N Gowda, Antonio H Ribeiro, Patrick Schwab, Kim Branson, Lei Clifton, Antonio Luiz P Ribeiro, Zhang-daihong Liu, and David A Clifton. Cardiac health assessment across scenarios and devices using a multimodal foundation model pretrained on data from 1.7 million individuals. *Nature Machine Intelligence*, 8:220–233, 2026. [1](#), [2](#), [13](#)
- [13] Yu Han, Vittorio Murino, Xiaofeng Liu, Xiang Zhang, and Cheng Ding. A systematic review on foundation models for electrocardiogram analysis: initial strides and expansive horizons. *arXiv preprint arXiv:2410.19877*, 2024. [1](#)
- [14] Aman Kansal, Emma Chen, Boyang Tom Jin, Pranav Rajpurkar, and David A Kim. MC-MED, multimodal clinical monitoring in the emergency department. *Scientific Data*, 12(1):1094, 2025. [1](#), [8](#), [13](#)
- [15] Aman Kansal, Emma Chen, Tom Jin, Pranav Rajpurkar, and David Kim. Multimodal Clinical Monitoring in the Emergency Department (MC-MED). *PhysioNet*, 2025. Version 1.0.1. [1](#), [8](#), [13](#)

- [16] Jun Li, Aaron D Aguirre, Valdey Moura Junior, Jiarui Jin, Che Liu, Lanhai Zhong, Chenxi Sun, Gari Clifford, M Brandon Westover, and Shenda Hong. An electrocardiogram foundation model built on over 10 million recordings. *NEJM AI*, 2(7):AIoa2401033, 2025. [1](#), [2](#), [13](#)
- [17] Yuxuan Liang, Haomin Wen, Yuqi Nie, Yushan Jiang, Ming Jin, Dongjin Song, Shirui Pan, and Qingsong Wen. Foundation models for time series analysis: A tutorial and survey. In *Proceedings of the 30th ACM SIGKDD conference on knowledge discovery and data mining*, pages 6555–6565, 2024. [1](#)
- [18] Riccardo Lunelli, Angus Nicolson, Samuel Martin Pröll, Sebastian Johannes Reinstadler, Axel Bauer, and Clemens Dlaska. BenchECG and xECG: a benchmark and baseline for ECG foundation models. *arXiv preprint arXiv:2509.10151*, 2025. [1](#), [2](#), [13](#)
- [19] Dominique Makowski, Tam Pham, Zen J. Lau, Jan C. Brammer, François Lespinasse, Hung Pham, Christopher Schölzel, and Annabel S H Chen. NeuroKit2: A python toolbox for neurophysiological signal processing, 2020. [9](#), [13](#)
- [20] Leland McInnes, John Healy, and James Melville. UMAP: Uniform manifold approximation and projection for dimension reduction. *arXiv preprint arXiv:1802.03426*, 2018. [11](#)
- [21] Kaden McKeen, Sameer Masood, Augustin Toma, Barry Rubin, and Bo Wang. ECG-FM: An open electrocardiogram foundation model. *Jamia Open*, 8(5):ooaf122, 2025. [1](#)
- [22] Fabian Pedregosa, Gaël Varoquaux, Alexandre Gramfort, Vincent Michel, Bertrand Thirion, Olivier Grisel, Mathieu Blondel, Peter Prettenhofer, Ron Weiss, Vincent Dubourg, Jake Vanderplas, Alexandre Passos, David Cournapeau, Matthieu Brucher, Matthieu Perrot, and Édouard Duchesnay. Scikit-learn: Machine learning in python. *Journal of Machine Learning Research (JMLR)*, 12(85):2825–2830, 2011. [12](#)
- [23] Manh Pham Hung, Aaqib Saeed, and Dong Ma. Boosting masked ECG-text auto-encoders as discriminative learners. In *International Conference on Machine Learning (ICML)*, pages 49277–49291, 2025. [1](#), [2](#), [13](#)
- [24] Arvind Pillai, Dimitris Spathis, Fahim Kawsar, and Mohammad Malekzadeh. PaPaGei: Open foundation models for optical physiological signals. In *International Conference on Learning Representations (ICLR)*, pages 48230–48261, 2025. [2](#), [13](#)
- [25] Mithun Saha, Maxwell A Xu, Wanting Mao, Sameer Neupane, James M Rehg, and Santosh Kumar. Pulse-PPG: An open-source field-trained PPG foundation model for wearable applications across lab and field settings. *Proceedings of the ACM on Interactive, Mobile, Wearable and Ubiquitous Technologies*, 9(3):1–35, 2025. [1](#)
- [26] Mohamed Adel Serhani, Hadeel T. El Kassabi, Heba Ismail, and Alramzana Nujum Navaz. ECG monitoring systems: Review, architecture, processes, and key challenges. *Sensors*, 20(6):1796, 2020. [1](#)
- [27] Sameer Sundrani, Julie Chen, Boyang Tom Jin, Zahra Shakeri Hossein Abad, Pranav Rajpurkar, and David Kim. Predicting patient decompensation from continuous physiologic monitoring in the emergency department. *npj Digital Medicine*, 6(1):60, 2023. [1](#)
- [28] Megha Thukral, Cyrus Tanade, Simon A Lee, Juhyeon Lee, Hao Zhou, Keum San Chun, Migyeong Gwak, Viswam Nathan, Md Mahbubur Rahman, Li Zhu, et al. Wavelet-driven masked multiscale reconstruction for PPG foundation models. *arXiv preprint arXiv:2601.12215*, 2026. [1](#)
- [29] Jiexia Ye, Yongzi Yu, Weiqi Zhang, Le Wang, Jia Li, and Fugee Tsung. Empowering time series analysis with foundation models: A comprehensive survey. *arXiv preprint arXiv:2405.02358*, 2024. [1](#)

**SignalMC-MED: A Multimodal Benchmark for Evaluating Biosignal
Foundation Models on Single-Lead ECG and PPG**

Supplementary Material

A. Supplementary Figures

This section contains Figure **S1 - S11**.

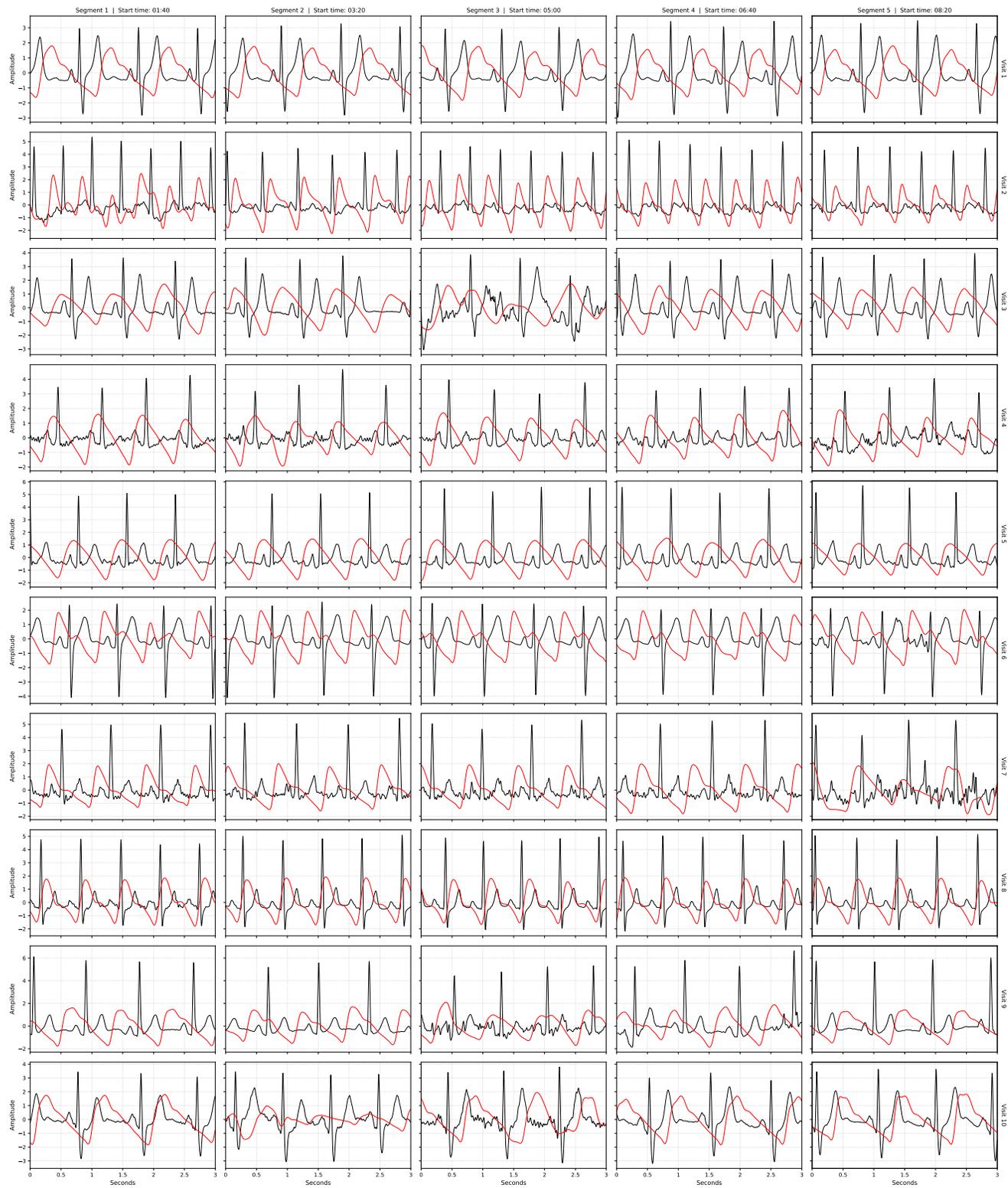


Figure S1. **Additional synchronized ECG and PPG signals in SignalMC-MED.** An extension of Figure 1b, with seven more examples of synchronized single-lead ECG (black) and PPG (red) signals from visits in the *train* set.

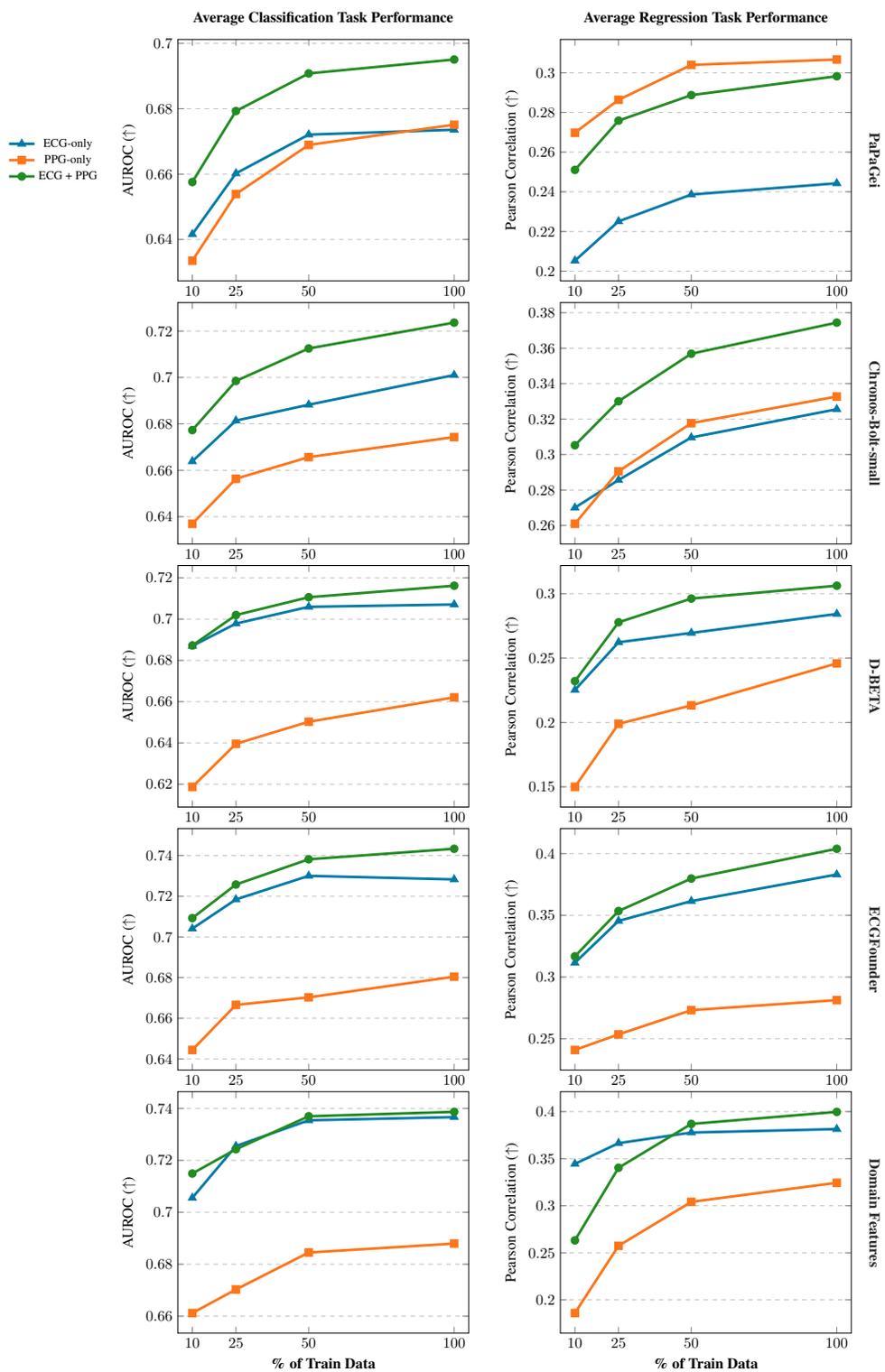


Figure S2. Multimodal fusion comparison on *val* for the remaining models. An extension of Figure 3, comparing the ECG-only, PPG-only, and ECG + PPG performance for the remaining models (PaPaGei, Chronos-Bolt-small, D-BETA, ECGFounder, Domain Features).

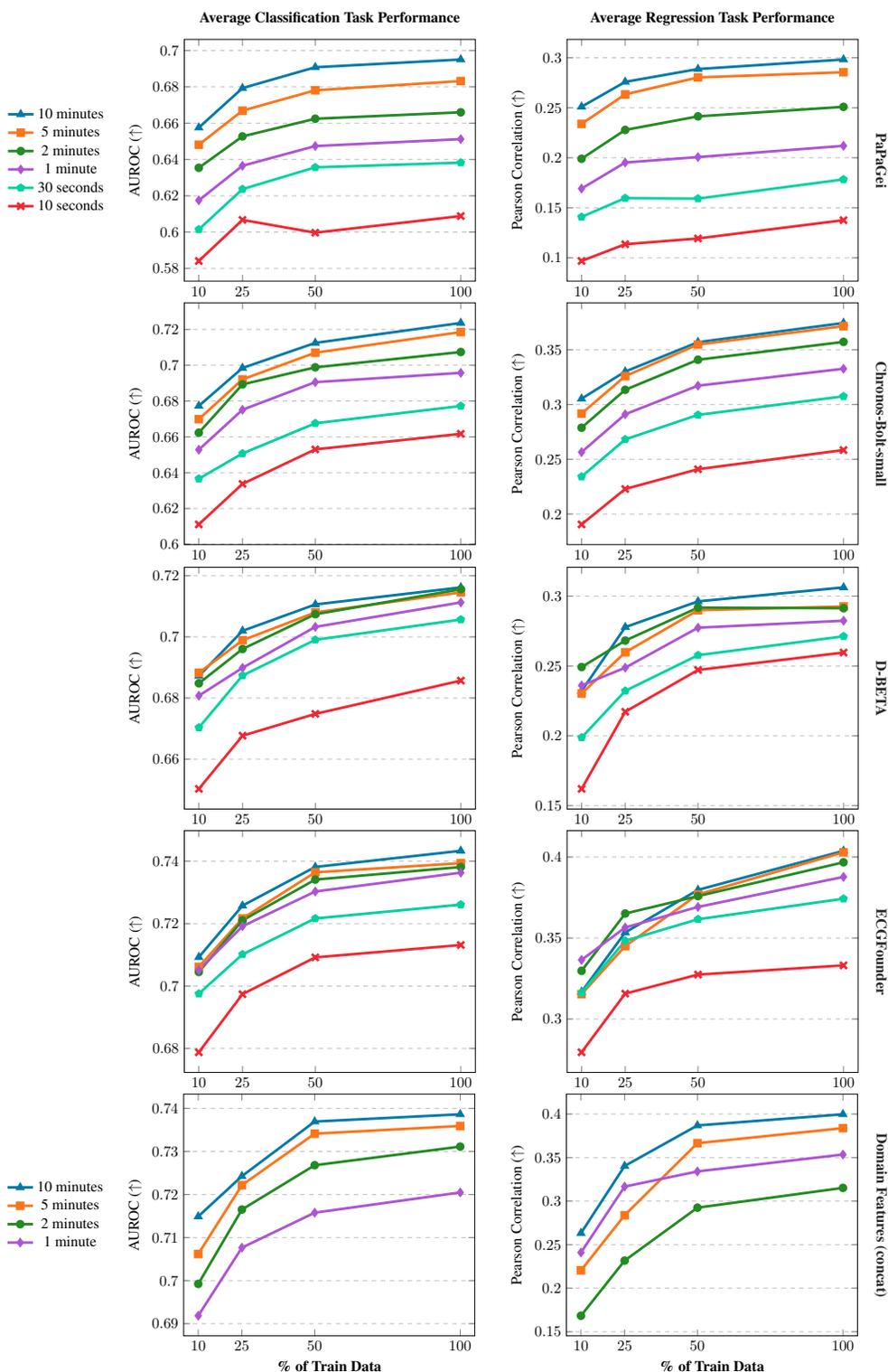


Figure S3. **Signal length comparison on val for the remaining models.** An extension of Figure 4, showing how the *ECG + PPG* performance is affected by increasing the signal length from 10 seconds to 10 minutes for the remaining models (PaPaGei, Chronos-Bolt-small, D-BETA, ECGFounder, Domain Features).

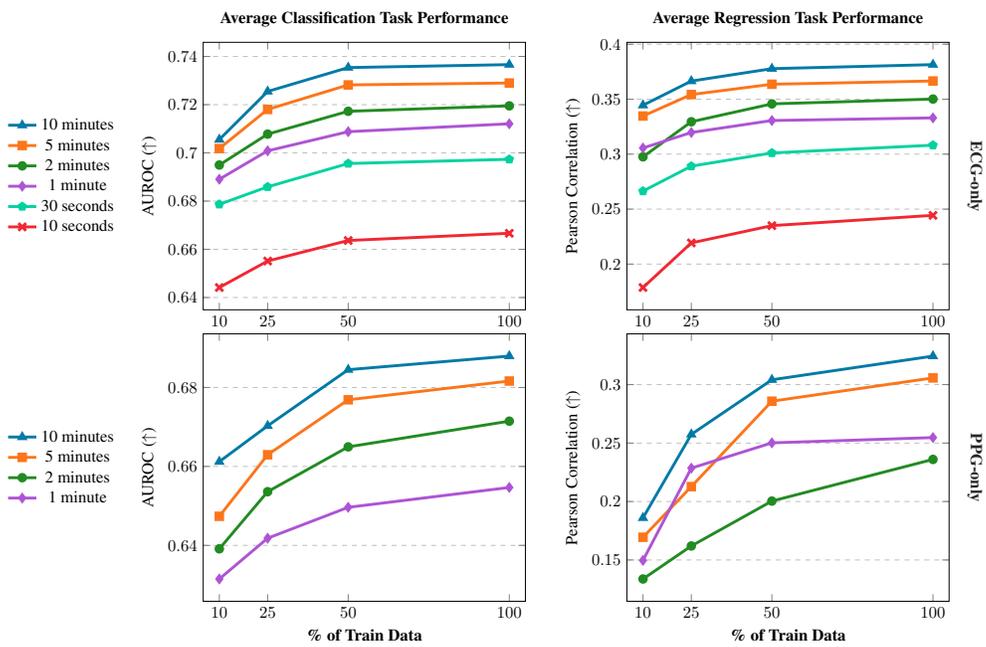


Figure S4. **Signal length comparison on val for hand-crafted ECG and PPG domain features**, showing how the performance is affected by increasing the signal length from 10 seconds to 10 minutes for ECG, and from 1 minute to 10 minutes for PPG. This analysis complements Figure 4 & S3.

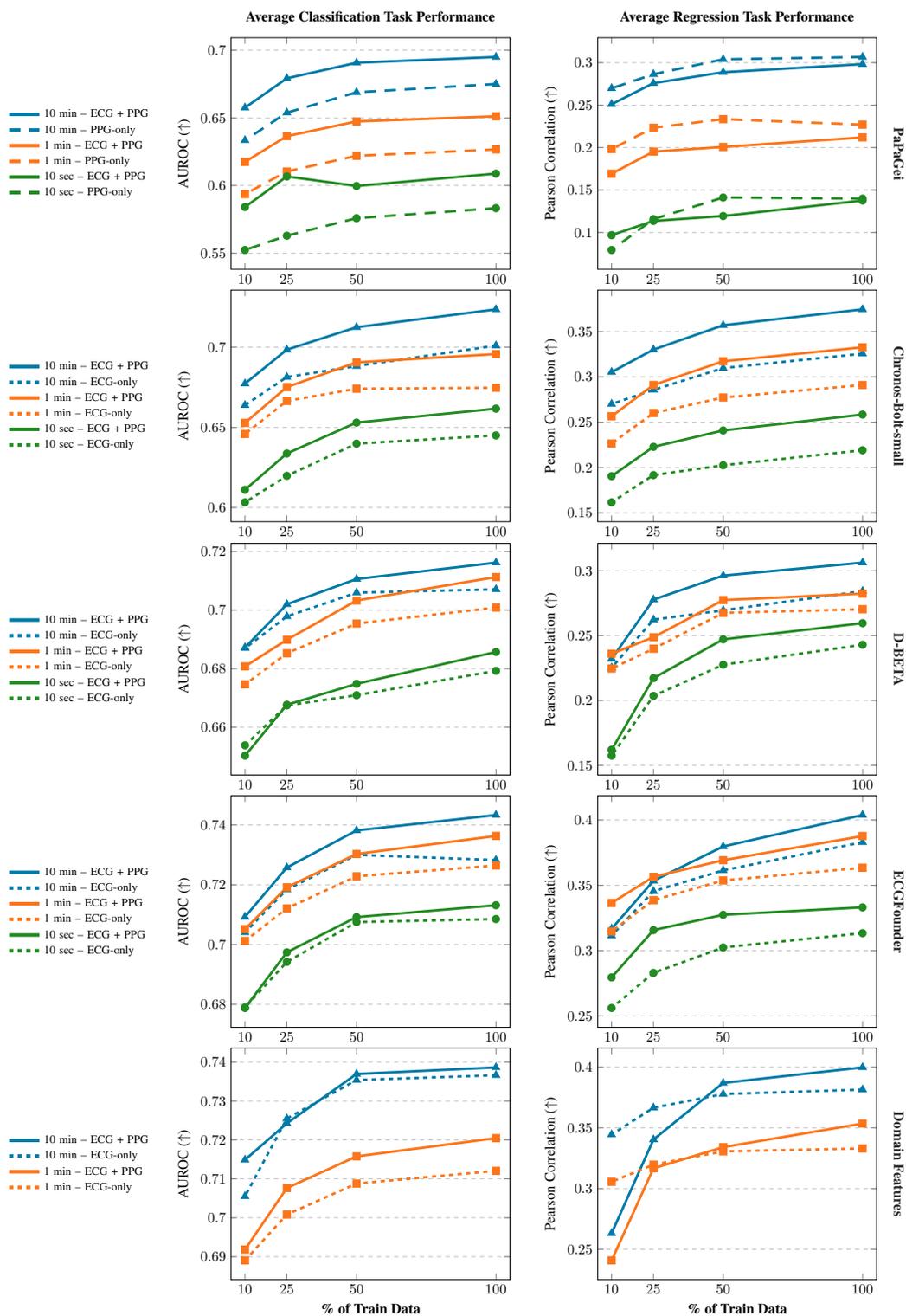


Figure S5. Multimodal fusion comparison across signal lengths on *val* for the remaining models. An extension of Figure 5, comparing the *PPG-only* and ECG + PPG performance for PaPaGei, and comparing the *ECG-only* and ECG + PPG performance for Chronos-Bolt-small, D-BETA, ECGFounder, and Domain Features.

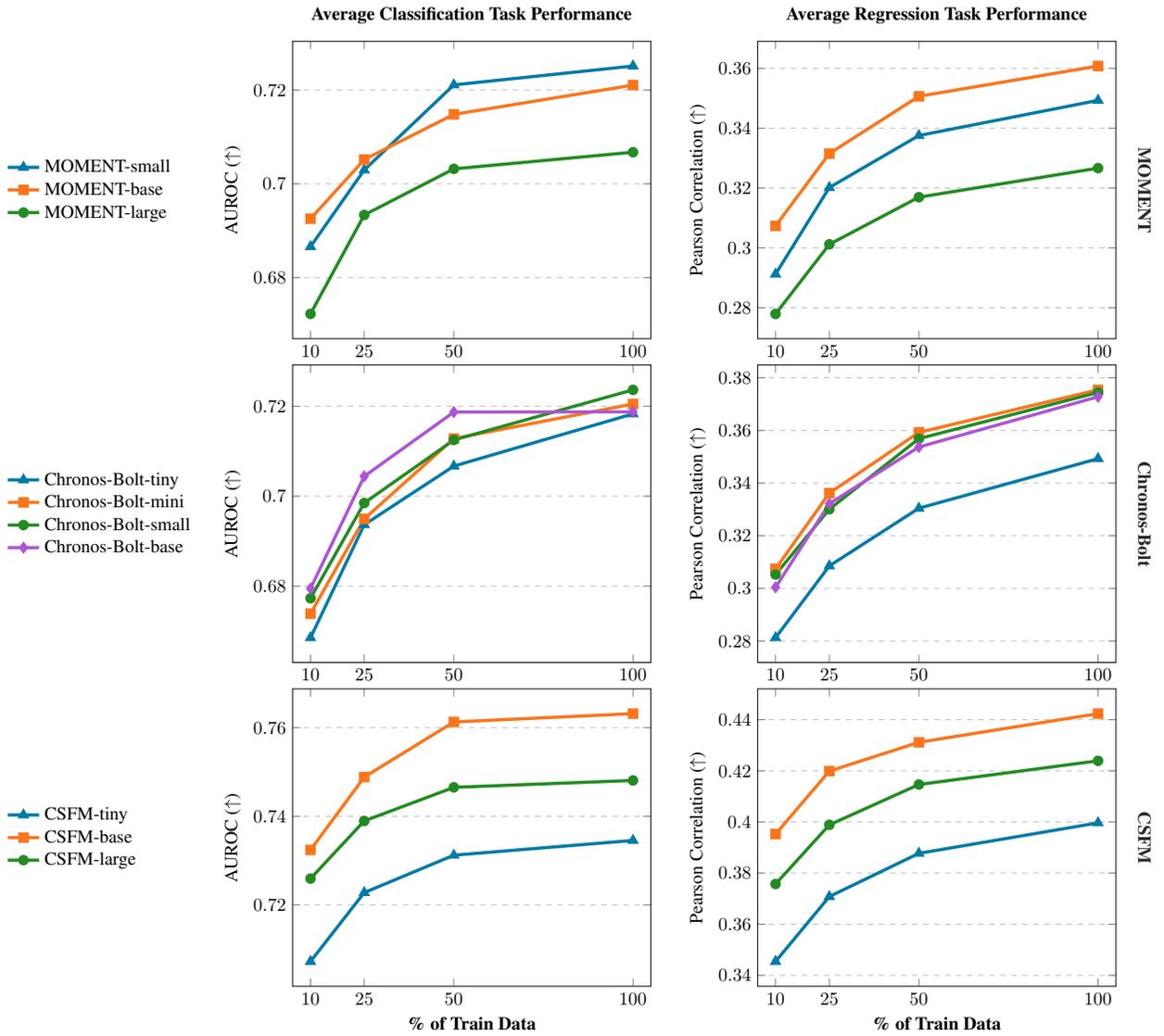


Figure S6. **Model size comparison on the val set**, comparing the *ECG + PPG* performance of small/base/large model variants for MOMENT, of tiny/mini/small/base variants for Chronos-Bolt, and of tiny/base/large variants for CSFM.

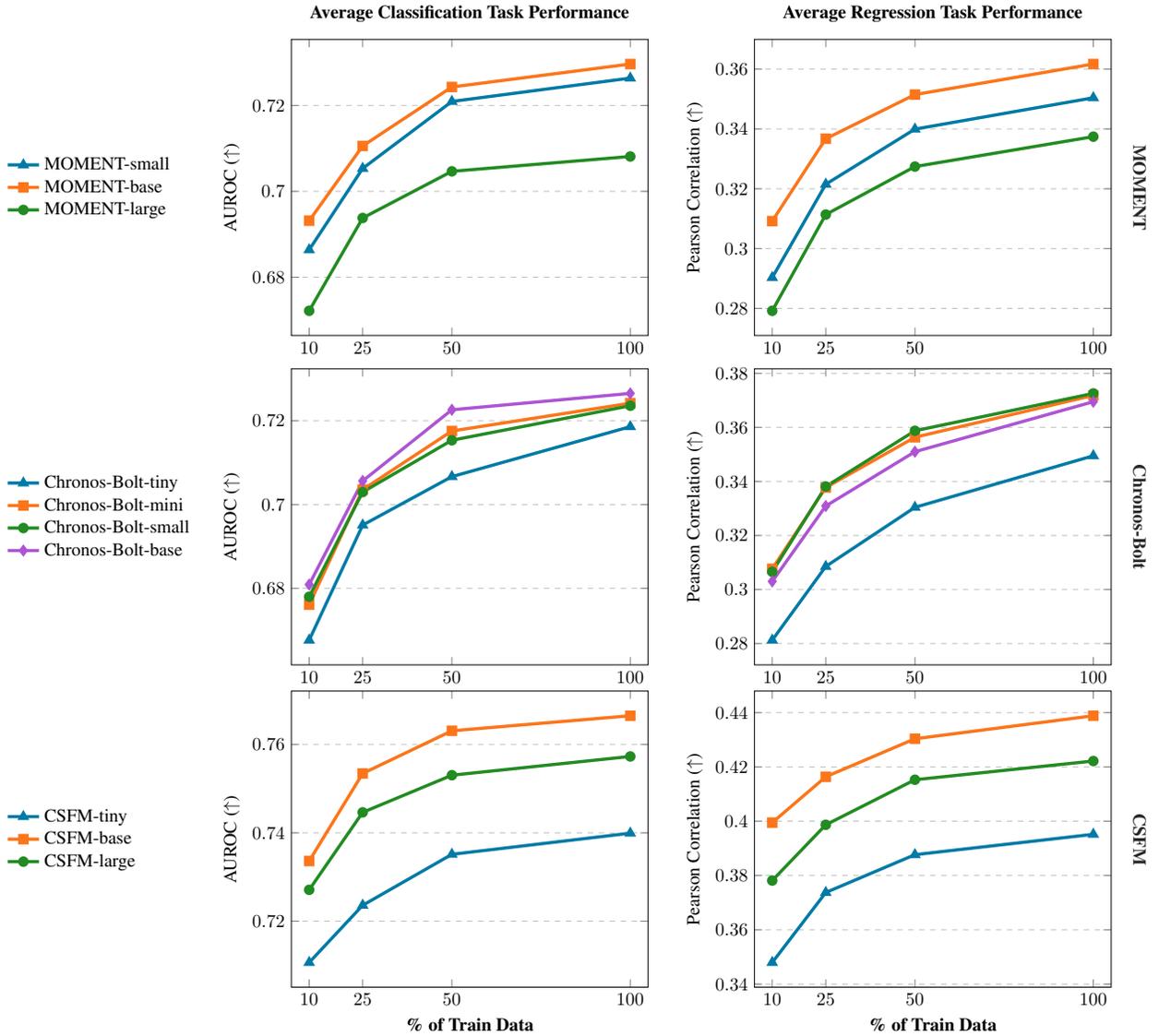


Figure S7. **Model size comparison with PCA on the *val* set.** The same analysis as Figure S6, but with all visit-level feature vectors projected down to 256 dimensions using PCA prior to downstream linear probing. Performance trends remain similar: larger model variants do *not* reliably outperform smaller ones. This indicates that differences in feature dimensionality are not the primary explanation for the occasional underperformance of larger models.

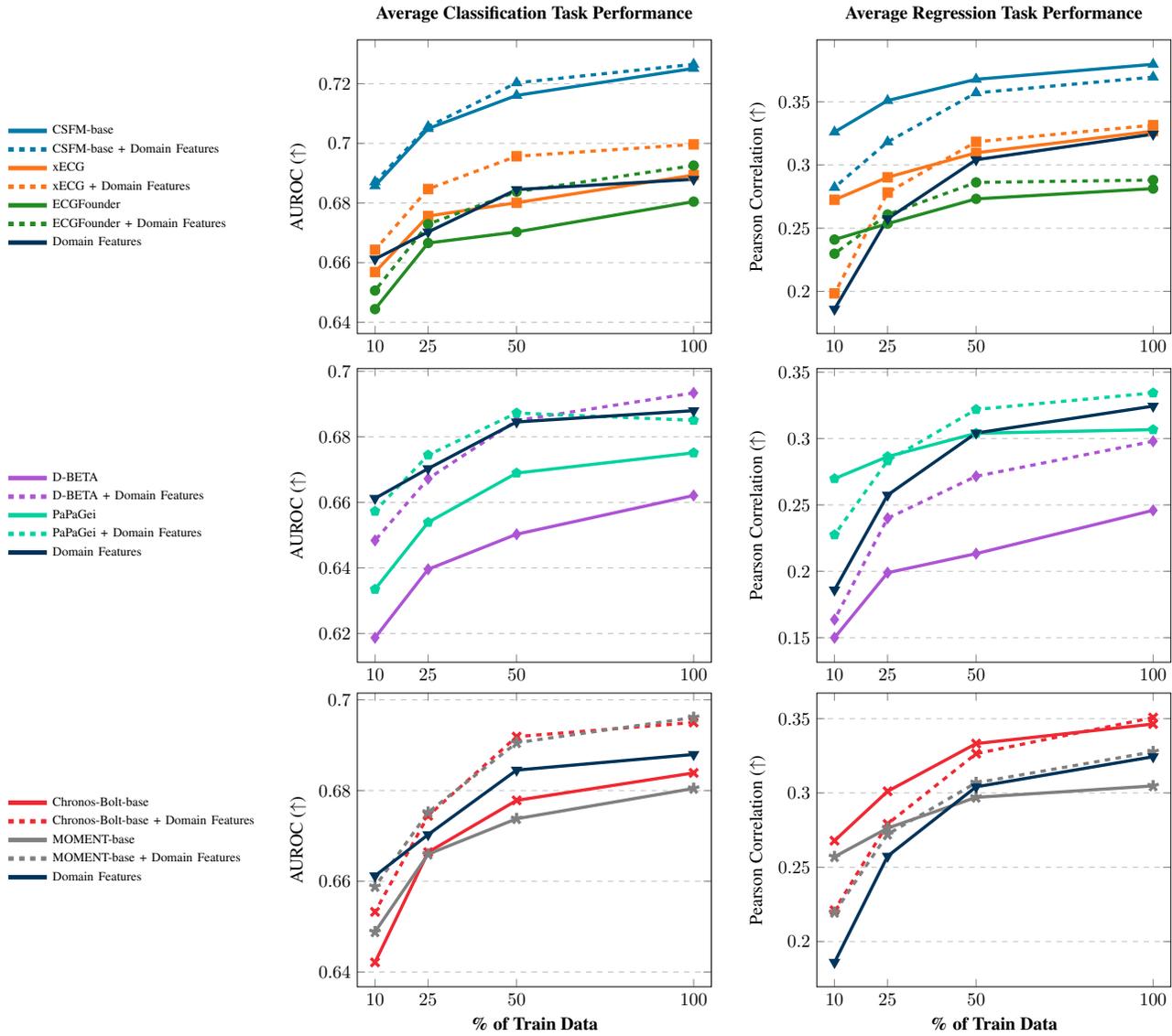


Figure S8. **Effect of concatenating PPG domain features with learned FM representations on the *val* set.** Comparison of PPG-only performance before (solid lines) and after (dashed lines) concatenation of hand-crafted PPG domain features to the frozen FM feature vectors, shown for each of the seven evaluated FMs. Identical analysis as in Figure 6, but for PPG domain features.

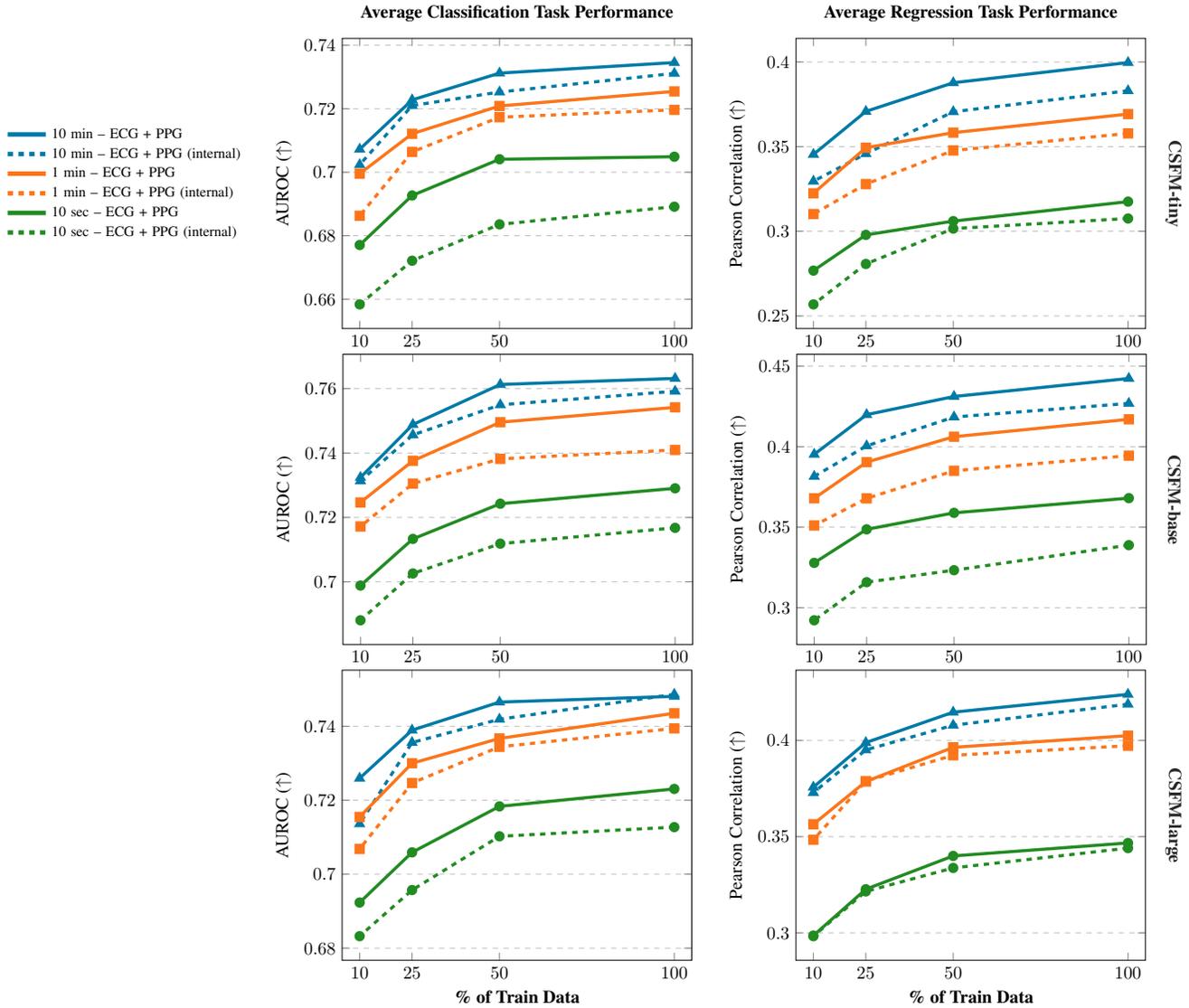


Figure S9. **Multimodal fusion comparison for CSFM on the *val* set.** Comparison between feature-level fusion (solid lines), obtained by averaging independently extracted ECG-only and PPG-only features, and CSFM’s internal multimodal fusion mechanism (dashed lines), which jointly processes ECG + PPG as a two-channel input. Results are shown for multiple signal lengths (10 seconds, 1 minute, 10 minutes) and all three CSFM model variants (tiny, base, large). Feature-level fusion consistently achieves slightly better performance.

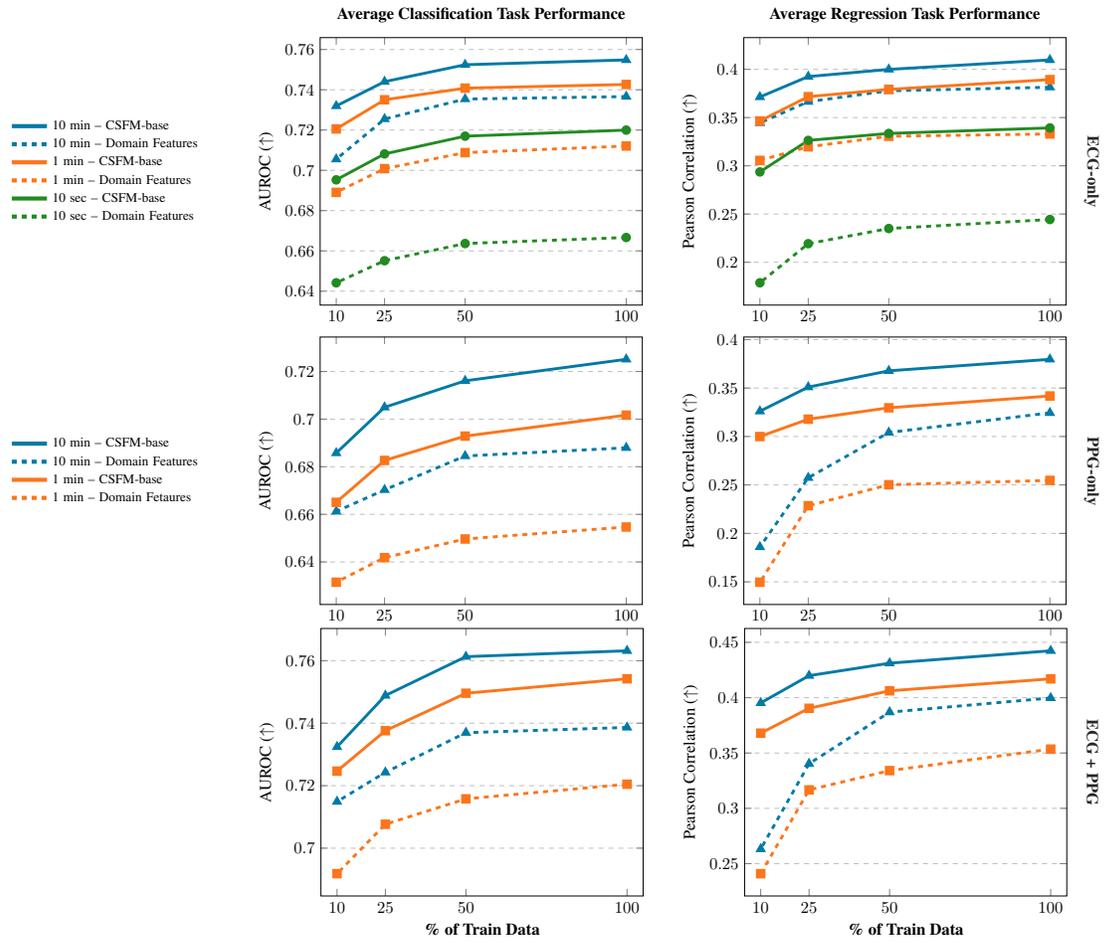


Figure S10. Comparison of CSFM-base and Domain Features across signal lengths on the *val* set, showing ECG-only, PPG-only, and ECG + PPG performance of CSFM-base (solid lines) and Domain Features (dashed lines) with multiple signal lengths. Especially for ECG-only, the performance gap narrows with longer signals.

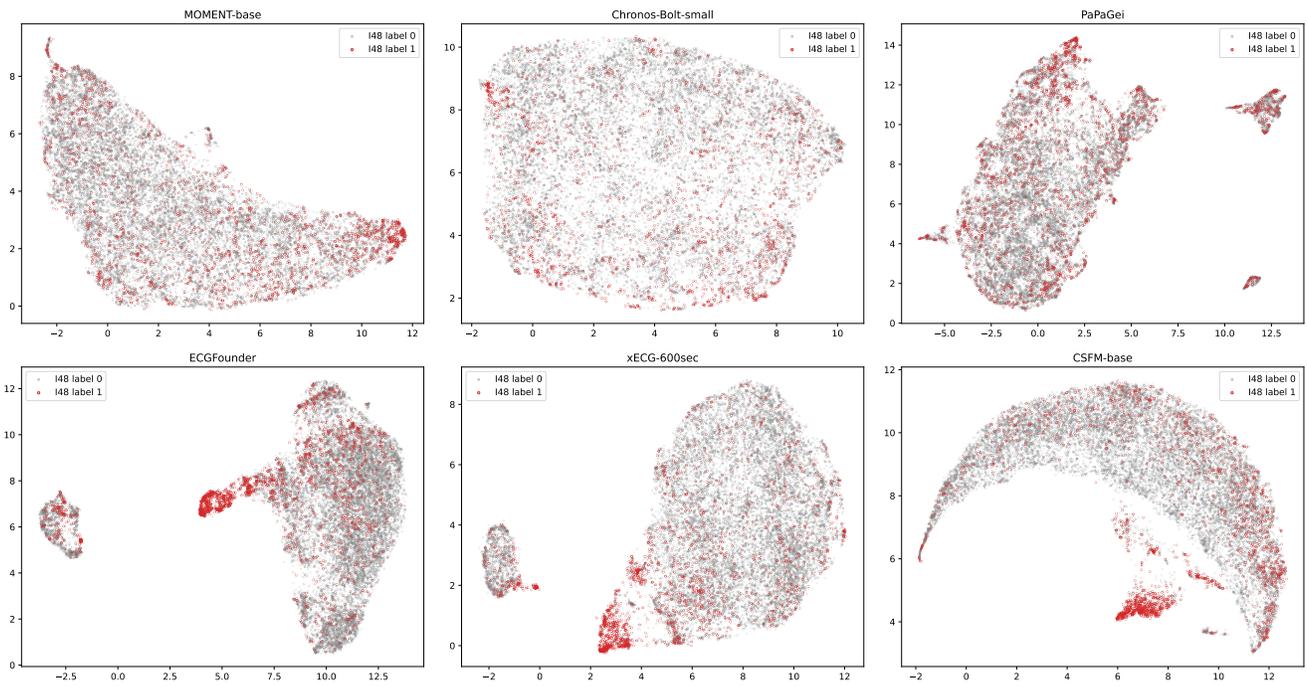


Figure S11. **UMAP visualizations colored by atrial fibrillation (I48) label.** UMAP projections of the *train* set *ECG + PPG* feature vectors for six selected models. CSFM-base, xECG-10min and ECGFounder exhibit clearer separation between positive and negative cases compared to general time-series FMs (MOMENT-base, Chronos-Bolt-small) and PaPaGei.

B. Supplementary Tables

This section contains Table **S1 - S16**.

Table S1. **Comparison of all model variants on the val set**, with results averaged across training with 10%, 25%, 50%, and 100% of the train visits (five resampling repeats per train fraction). Regression tasks are evaluated using Pearson correlation (\uparrow), and classification tasks using AUROC (\uparrow). **Bold red** marks the best result in each column, **bold** marks the second-best, and underline marks the third-best. Based on these results and the corresponding ranking in Table S2, one variant each for MOMENT, Chronos-Bolt, xECG, CSFM, and the PPG domain features is selected for inclusion in the main model comparison on the test set. These selected variants are here indicated with italics and a superscript * next to the model name.

Model Name	Modality	Age reg	Sex cls	ED.dispo cls	Labs reg	ICD-10 cls
MOMENT-small	ECG-only	0.668	0.750	0.696	0.239	0.685
MOMENT-small	PPG-only	0.676	0.721	0.665	0.224	0.652
MOMENT-small	ECG + PPG	0.747	0.781	0.706	0.272	0.701
<i>MOMENT-base*</i>	ECG-only	0.680	0.757	0.695	0.253	0.683
<i>MOMENT-base*</i>	PPG-only	0.691	0.725	0.680	0.233	0.660
<i>MOMENT-base*</i>	ECG + PPG	0.762	0.791	0.703	0.284	0.700
MOMENT-large	ECG-only	0.635	0.720	0.687	0.223	0.673
MOMENT-large	PPG-only	0.638	0.690	0.645	0.200	0.628
MOMENT-large	ECG + PPG	0.722	0.755	0.690	0.254	0.688
Chronos-Bolt-tiny	ECG-only	0.667	0.721	0.687	0.236	0.685
Chronos-Bolt-tiny	PPG-only	0.680	0.718	0.682	0.233	0.649
Chronos-Bolt-tiny	ECG + PPG	0.739	0.761	0.696	0.265	0.690
Chronos-Bolt-mini	ECG-only	0.677	0.732	0.697	0.252	0.671
Chronos-Bolt-mini	PPG-only	0.707	0.721	0.681	0.248	0.660
Chronos-Bolt-mini	ECG + PPG	0.762	0.777	0.707	0.292	0.691
<i>Chronos-Bolt-small*</i>	ECG-only	0.687	0.745	0.696	0.249	0.675
<i>Chronos-Bolt-small*</i>	PPG-only	0.705	0.719	0.676	0.250	0.650
<i>Chronos-Bolt-small*</i>	ECG + PPG	0.765	0.785	0.707	0.289	0.693
Chronos-Bolt-base	ECG-only	0.670	0.749	0.694	0.245	0.685
Chronos-Bolt-base	PPG-only	0.726	0.725	0.688	0.260	0.659
Chronos-Bolt-base	ECG + PPG	0.764	0.784	0.704	0.287	0.697
xECG	ECG-only	0.768	0.805	0.728	0.323	0.721
xECG	PPG-only	0.737	0.750	0.685	0.245	0.666
xECG	ECG + PPG	0.815	<u>0.829</u>	0.731	0.346	0.726
<i>xECG-10min*</i>	ECG-only	0.771	0.809	0.730	0.334	0.724
<i>xECG-10min*</i>	PPG-only	0.756	0.755	0.695	0.264	0.672
<i>xECG-10min*</i>	ECG + PPG	<u>0.818</u>	0.833	<u>0.733</u>	0.352	0.729
ECGFounder	ECG-only	0.760	0.769	0.716	0.299	0.715
ECGFounder	PPG-only	0.672	0.724	0.674	0.211	0.658
ECGFounder	ECG + PPG	0.808	0.804	0.721	0.308	0.722
D-BETA	ECG-only	0.631	0.719	0.692	0.214	0.698
D-BETA	PPG-only	0.599	0.700	0.657	0.152	0.635
D-BETA	ECG + PPG	0.675	0.755	0.695	0.228	0.699
PaPaGei	ECG-only	0.542	0.670	0.648	0.189	0.663
PaPaGei	PPG-only	0.664	0.688	0.656	0.245	0.655
PaPaGei	ECG + PPG	0.676	0.722	0.665	0.229	0.678
CSFM-tiny	ECG-only	0.703	0.774	0.705	0.278	0.706
CSFM-tiny	PPG-only	0.746	0.750	0.693	0.270	0.671
CSFM-tiny	ECG + PPG	0.800	0.809	0.721	0.323	0.715
<i>CSFM-base*</i>	ECG-only	0.775	0.804	0.735	0.346	0.741
<i>CSFM-base*</i>	PPG-only	0.784	0.758	0.720	0.303	0.701
<i>CSFM-base*</i>	ECG + PPG	0.834	<u>0.829</u>	0.740	0.371	0.744
CSFM-large	ECG-only	0.761	0.810	0.721	0.320	0.720
CSFM-large	PPG-only	0.759	0.763	0.700	0.268	0.686
CSFM-large	ECG + PPG	0.824	0.839	0.732	<u>0.351</u>	0.730
ECG Domain Features	ECG-only	0.690	0.668	0.723	0.327	<u>0.732</u>
PPG Domain Features - 20sec	PPG-only	0.643	0.731	0.693	0.190	0.668
<i>PPG Domain Features - 60sec*</i>	PPG-only	0.699	0.732	0.694	0.214	0.668
PPG Domain Features - 120sec	PPG-only	0.690	0.733	0.695	0.204	0.669
ECG Domain Features + PPG Domain Features - 20sec (concat)	ECG + PPG	0.745	0.760	0.728	0.271	<u>0.732</u>
<i>ECG Domain Features + PPG Domain Features - 60sec (concat)*</i>	ECG + PPG	0.784	0.760	0.725	0.293	0.726
ECG Domain Features + PPG Domain Features - 120sec (concat)	ECG + PPG	0.773	0.761	0.726	0.275	0.729

Table S2. **Ranking of all model variants on the *val* set**, computed jointly across ECG-only, PPG-only, and ECG + PPG. Mean model rank (\downarrow) across the five aggregated task categories based on Table S1. Based on this ranking, one variant each for MOMENT, Chronos-Bolt, xECG, CSFM, and the PPG domain features is selected for inclusion in the main model comparison on the *test* set.

Rank	Model Name	Modality	Model Ranks	Mean Model Rank (\downarrow)
1	CSFM-base	ECG + PPG	1,4,1,1,1	1.6
2	CSFM-large	ECG + PPG	2,1,4,3,5	3.0
3	xECG-10min	ECG + PPG	3,2,3,2,7	3.4
4	xECG	ECG + PPG	4,3,5,4,8	4.8
5	CSFM-base	ECG-only	9,9,2,5,2	5.4
6	xECG-10min	ECG-only	11,6,6,6,10	7.8
7	xECG	ECG-only	12,8,7,8,12	9.4
8	CSFM-tiny	ECG + PPG	6,7,12,9,15	9.8
9	ECGFounder	ECG + PPG	5,10,13,11,11	10.0
10	CSFM-large	ECG-only	17,5,14,10,13	11.8
11	ECG Domain Features + PPG Domain Features - 60sec (concat)	ECG + PPG	8,21,10,14,9	12.4
12	ECG Domain Features + PPG Domain Features - 120sec (concat)	ECG + PPG	10,19,9,20,6	12.8
13	CSFM-base	PPG-only	7,23,15,12,18	15.0
14	ECGFounder	ECG-only	18,17,16,13,14	15.6
15	ECG Domain Features + PPG Domain Features - 20sec (concat)	ECG + PPG	23,22,8,22,4	15.8
16	Chronos-Bolt-small	ECG + PPG	13,12,18,16,23	16.4
17	MOMENT-base	ECG + PPG	15,11,22,18,19	17.0
18	Chronos-Bolt-mini	ECG + PPG	16,15,17,15,24	17.4
19	Chronos-Bolt-base	ECG + PPG	14,13,21,17,22	17.4
20	MOMENT-small	ECG + PPG	21,14,19,21,17	18.4
21	CSFM-tiny	ECG-only	30,16,20,19,16	20.2
22	ECG Domain Features	ECG-only	33,52,11,7,3	21.2
23	CSFM-large	PPG-only	19,18,23,24,27	22.2
24	Chronos-Bolt-tiny	ECG + PPG	24,20,27,25,25	24.2
25	xECG-10min	PPG-only	20,26,30,26,35	27.4
26	MOMENT-large	ECG + PPG	27,25,37,28,26	28.6
27	CSFM-tiny	PPG-only	22,29,35,23,36	29.0
28	MOMENT-base	ECG-only	37,24,28,29,31	29.8
29	Chronos-Bolt-small	ECG-only	35,32,26,32,33	31.6
30	D-BETA	ECG + PPG	41,27,29,42,20	31.8
31	MOMENT-small	ECG-only	44,28,25,37,30	32.8
32	Chronos-Bolt-mini	ECG-only	38,35,24,30,37	32.8
33	Chronos-Bolt-base	ECG-only	43,31,33,36,28	34.2
34	xECG	PPG-only	25,30,41,35,41	34.4
35	Chronos-Bolt-base	PPG-only	26,37,38,27,45	34.6
36	PPG Domain Features - 60sec	PPG-only	31,34,32,45,39	36.2
37	PPG Domain Features - 120sec	PPG-only	34,33,31,48,38	36.8
38	Chronos-Bolt-mini	PPG-only	28,43,43,33,43	38.0
39	Chronos-Bolt-tiny	ECG-only	45,41,39,38,29	38.4
40	MOMENT-base	PPG-only	32,38,44,39,44	39.4
41	D-BETA	ECG-only	50,45,36,46,21	39.6
42	Chronos-Bolt-small	PPG-only	29,46,45,31,49	40.0
43	PaPaGei	ECG + PPG	40,40,47,41,32	40.0
44	PPG Domain Features - 20sec	PPG-only	47,36,34,50,40	41.4
45	MOMENT-large	ECG-only	49,44,40,44,34	42.2
46	Chronos-Bolt-tiny	PPG-only	36,47,42,40,50	43.0
47	MOMENT-small	PPG-only	39,42,48,43,48	44.0
48	ECGFounder	PPG-only	42,39,46,47,46	44.0
49	PaPaGei	PPG-only	46,50,50,34,47	45.4
50	PaPaGei	ECG-only	52,51,51,51,42	49.4
51	MOMENT-large	PPG-only	48,49,52,49,52	50.0
52	D-BETA	PPG-only	51,48,49,52,51	50.2

Table S3. **Main model ranking on the *test* set**, computed jointly across ECG-only, PPG-only, and ECG + PPG. Mean model rank (\downarrow) across the five aggregated task categories based on Table 2. This table complements the modality-specific rankings shown in Table 3.

Rank	Model Name	Modality	Model Ranks	Mean Model Rank (\downarrow)
1	CSFM-base	ECG + PPG	1,1,1,1,3	1.4
2	xECG-10min	ECG + PPG	2,2,6,3,6	3.8
3	CSFM-base	ECG-only	7,5,5,2,4	4.6
4	ECGFounder	ECG + PPG	3,4,8,5,5	5.0
5	Domain Features (concat)	ECG + PPG	8,7,2,9,2	5.6
6	xECG-10min	ECG-only	10,3,7,4,8	6.4
7	CSFM-base	PPG-only	4,9,3,8,9	6.6
8	MOMENT-base	ECG + PPG	5,6,11,10,11	8.6
9	ECGFounder	ECG-only	11,12,9,7,7	9.2
10	Domain Features	ECG-only	17,23,4,6,1	10.2
11	xECG-10min	PPG-only	6,10,10,12,14	10.4
12	Chronos-Bolt-small	ECG + PPG	9,8,16,11,13	11.4
13	D-BETA	ECG + PPG	14,14,15,14,10	13.4
14	MOMENT-base	PPG-only	12,17,12,13,17	14.2
15	Domain Features	PPG-only	15,13,13,22,15	15.6
16	MOMENT-base	ECG-only	21,11,19,15,16	16.4
17	Chronos-Bolt-small	PPG-only	13,15,17,16,22	16.6
18	ECGFounder	PPG-only	16,16,14,19,18	16.6
19	D-BETA	ECG-only	22,21,18,20,12	18.6
20	Chronos-Bolt-small	ECG-only	19,18,20,18,19	18.8
21	PaPaGei	PPG-only	18,22,22,17,21	20.0
22	PaPaGei	ECG + PPG	20,19,23,21,20	20.6
23	D-BETA	PPG-only	23,20,21,24,23	22.2
24	PaPaGei	ECG-only	24,24,24,23,24	23.8

Table S4. **Model comparison of ECG-only *val* performance on all 20 benchmark tasks when using 100% of the train visits**. Tasks are grouped into two demographics tasks, ED disposition classification, eight laboratory value regression tasks, and nine prior ICD-10 diagnosis classification tasks. Performance on the *val* set when using 100% of the train set visits. Training visit are sampled with replacement and both downstream training and hyperparameter selection are repeated five times, results are mean \pm std across repetitions. Regression tasks are evaluated using Pearson correlation (\uparrow), and classification tasks using AUROC (\uparrow). **Red red** marks the best mean value in each row.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.696 \pm 0.001	0.694 \pm 0.001	0.710 \pm 0.001	0.559 \pm 0.001	0.688 \pm 0.005	0.778 \pm 0.003	0.785 \pm 0.001	0.785 \pm 0.002
Sex cls	0.676 \pm 0.001	0.774 \pm 0.002	0.770 \pm 0.003	0.675 \pm 0.002	0.730 \pm 0.003	0.788 \pm 0.003	0.822 \pm 0.002	0.816 \pm 0.002
ED_dispo cls	0.729 \pm 0.001	0.703 \pm 0.003	0.712 \pm 0.005	0.655 \pm 0.003	0.695 \pm 0.004	0.725 \pm 0.003	0.745 \pm 0.001	0.744 \pm 0.001
Potassium reg	0.268 \pm 0.004	0.322 \pm 0.008	0.304 \pm 0.010	0.194 \pm 0.016	0.188 \pm 0.023	0.315 \pm 0.005	0.383 \pm 0.007	0.386 \pm 0.007
Calcium reg	0.155 \pm 0.011	0.150 \pm 0.007	0.127 \pm 0.018	0.077 \pm 0.012	0.098 \pm 0.016	0.155 \pm 0.014	0.205 \pm 0.010	0.209 \pm 0.010
eGFR reg	0.546 \pm 0.001	0.494 \pm 0.007	0.509 \pm 0.005	0.369 \pm 0.006	0.481 \pm 0.015	0.579 \pm 0.007	0.585 \pm 0.003	0.591 \pm 0.006
Glucose reg	0.297 \pm 0.007	0.177 \pm 0.010	0.199 \pm 0.008	0.157 \pm 0.013	0.150 \pm 0.018	0.244 \pm 0.009	0.257 \pm 0.014	0.262 \pm 0.017
Hemoglobin reg	0.359 \pm 0.002	0.267 \pm 0.007	0.261 \pm 0.008	0.205 \pm 0.014	0.253 \pm 0.014	0.356 \pm 0.015	0.365 \pm 0.008	0.373 \pm 0.007
Albumin reg	0.434 \pm 0.003	0.290 \pm 0.006	0.318 \pm 0.007	0.278 \pm 0.006	0.262 \pm 0.023	0.414 \pm 0.006	0.422 \pm 0.006	0.428 \pm 0.004
BUN reg	0.388 \pm 0.006	0.322 \pm 0.004	0.332 \pm 0.014	0.224 \pm 0.015	0.289 \pm 0.010	0.376 \pm 0.013	0.384 \pm 0.011	0.394 \pm 0.005
Sodium reg	0.291 \pm 0.009	0.189 \pm 0.015	0.171 \pm 0.010	0.134 \pm 0.010	0.150 \pm 0.014	0.228 \pm 0.013	0.249 \pm 0.010	0.260 \pm 0.009
I48 cls	0.826 \pm 0.003	0.821 \pm 0.009	0.805 \pm 0.004	0.761 \pm 0.006	0.838 \pm 0.007	0.844 \pm 0.009	0.845 \pm 0.009	0.861 \pm 0.005
I50 cls	0.846 \pm 0.002	0.812 \pm 0.010	0.812 \pm 0.018	0.769 \pm 0.003	0.826 \pm 0.005	0.819 \pm 0.008	0.842 \pm 0.002	0.840 \pm 0.004
{Z95.0, Z95.2, Z95.3} cls	0.831 \pm 0.003	0.778 \pm 0.018	0.789 \pm 0.017	0.736 \pm 0.006	0.819 \pm 0.020	0.833 \pm 0.012	0.825 \pm 0.009	0.854 \pm 0.010
N18 cls	0.798 \pm 0.003	0.707 \pm 0.019	0.727 \pm 0.014	0.709 \pm 0.009	0.730 \pm 0.003	0.749 \pm 0.004	0.775 \pm 0.004	0.783 \pm 0.005
{E10, E11} cls	0.707 \pm 0.004	0.628 \pm 0.013	0.625 \pm 0.008	0.623 \pm 0.010	0.640 \pm 0.005	0.678 \pm 0.009	0.695 \pm 0.005	0.696 \pm 0.005
{G47.3, E66.2} cls	0.677 \pm 0.010	0.649 \pm 0.022	0.635 \pm 0.018	0.633 \pm 0.015	0.648 \pm 0.013	0.671 \pm 0.017	0.696 \pm 0.011	0.686 \pm 0.016
{D50, D51, D52, D53, D64} cls	0.673 \pm 0.003	0.603 \pm 0.005	0.595 \pm 0.009	0.615 \pm 0.004	0.610 \pm 0.007	0.648 \pm 0.009	0.664 \pm 0.004	0.671 \pm 0.010
{J96, R09.02, R06.0} cls	0.678 \pm 0.002	0.636 \pm 0.008	0.636 \pm 0.008	0.630 \pm 0.006	0.654 \pm 0.007	0.660 \pm 0.010	0.680 \pm 0.005	0.676 \pm 0.004
{I26, I80, I82} cls	0.660 \pm 0.005	0.608 \pm 0.011	0.604 \pm 0.013	0.605 \pm 0.015	0.587 \pm 0.010	0.597 \pm 0.010	0.659 \pm 0.012	0.675 \pm 0.011

Table S5. Model comparison of PPG-only *val* performance on all 20 benchmark tasks when using 100% of the train visits. Identical format and content as Table S4, but for PPG-only. Performance on the *val* set when using 100% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.734 ± 0.003	0.704 ± 0.002	0.726 ± 0.002	0.673 ± 0.003	0.663 ± 0.007	0.703 ± 0.005	0.778 ± 0.002	0.799 ± 0.002
Sex cls	0.744 ± 0.004	0.731 ± 0.004	0.733 ± 0.005	0.691 ± 0.002	0.715 ± 0.004	0.735 ± 0.004	0.767 ± 0.003	0.767 ± 0.003
ED_dispo cls	0.706 ± 0.003	0.688 ± 0.005	0.689 ± 0.001	0.673 ± 0.005	0.671 ± 0.006	0.686 ± 0.003	0.707 ± 0.004	0.726 ± 0.001
Potassium reg	0.152 ± 0.016	0.175 ± 0.016	0.178 ± 0.008	0.164 ± 0.014	0.115 ± 0.025	0.135 ± 0.025	0.145 ± 0.007	0.190 ± 0.006
Calcium reg	0.149 ± 0.010	0.122 ± 0.009	0.156 ± 0.015	0.117 ± 0.012	0.100 ± 0.009	0.121 ± 0.040	0.165 ± 0.016	0.183 ± 0.015
eGFR reg	0.485 ± 0.005	0.439 ± 0.004	0.479 ± 0.004	0.436 ± 0.007	0.399 ± 0.014	0.419 ± 0.063	0.513 ± 0.007	0.567 ± 0.005
Glucose reg	0.281 ± 0.014	0.255 ± 0.004	0.292 ± 0.013	0.236 ± 0.012	0.193 ± 0.017	0.188 ± 0.062	0.266 ± 0.015	0.297 ± 0.009
Hemoglobin reg	0.321 ± 0.008	0.323 ± 0.008	0.311 ± 0.009	0.321 ± 0.012	0.222 ± 0.008	0.310 ± 0.024	0.338 ± 0.011	0.392 ± 0.009
Albumin reg	0.325 ± 0.002	0.289 ± 0.005	0.331 ± 0.010	0.317 ± 0.008	0.199 ± 0.022	0.278 ± 0.035	0.331 ± 0.011	0.394 ± 0.006
BUN reg	0.282 ± 0.008	0.235 ± 0.017	0.299 ± 0.011	0.282 ± 0.011	0.188 ± 0.020	0.223 ± 0.037	0.310 ± 0.022	0.368 ± 0.004
Sodium reg	0.191 ± 0.020	0.199 ± 0.006	0.224 ± 0.012	0.214 ± 0.007	0.136 ± 0.050	0.155 ± 0.028	0.225 ± 0.012	0.227 ± 0.012
I48 cls	0.767 ± 0.010	0.774 ± 0.002	0.758 ± 0.006	0.769 ± 0.002	0.753 ± 0.010	0.780 ± 0.003	0.780 ± 0.005	0.814 ± 0.004
I50 cls	0.731 ± 0.010	0.727 ± 0.004	0.732 ± 0.005	0.739 ± 0.006	0.689 ± 0.006	0.726 ± 0.006	0.747 ± 0.004	0.773 ± 0.006
{Z95.0, Z95.2, Z95.3} cls	0.656 ± 0.017	0.682 ± 0.011	0.671 ± 0.012	0.674 ± 0.008	0.674 ± 0.011	0.706 ± 0.014	0.709 ± 0.014	0.735 ± 0.016
N18 cls	0.698 ± 0.014	0.662 ± 0.005	0.676 ± 0.009	0.672 ± 0.007	0.676 ± 0.010	0.681 ± 0.015	0.714 ± 0.005	0.748 ± 0.004
{E10, E11} cls	0.684 ± 0.008	0.650 ± 0.006	0.653 ± 0.008	0.631 ± 0.008	0.651 ± 0.009	0.668 ± 0.007	0.674 ± 0.012	0.702 ± 0.012
{G47.3, E66.2} cls	0.662 ± 0.019	0.660 ± 0.005	0.664 ± 0.005	0.674 ± 0.010	0.629 ± 0.012	0.662 ± 0.007	0.684 ± 0.011	0.715 ± 0.006
{D50, D51, D52, D53, D64} cls	0.642 ± 0.006	0.620 ± 0.008	0.602 ± 0.008	0.615 ± 0.007	0.619 ± 0.008	0.606 ± 0.010	0.643 ± 0.014	0.665 ± 0.009
{J96, R09.02, R06.0} cls	0.667 ± 0.005	0.638 ± 0.004	0.634 ± 0.005	0.653 ± 0.007	0.624 ± 0.008	0.626 ± 0.004	0.638 ± 0.003	0.673 ± 0.004
{I26, I80, I82} cls	0.610 ± 0.017	0.654 ± 0.010	0.605 ± 0.018	0.636 ± 0.007	0.583 ± 0.008	0.608 ± 0.021	0.635 ± 0.013	0.659 ± 0.011

Table S6. Model comparison of ECG + PPG *val* performance on all 20 benchmark tasks when using 100% of the train visits. Identical format and content as Table S4, but for ECG + PPG. Performance on the *val* set when using 100% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.803 ± 0.003	0.776 ± 0.001	0.781 ± 0.002	0.697 ± 0.002	0.747 ± 0.003	0.822 ± 0.002	0.831 ± 0.001	0.842 ± 0.001
Sex cls	0.771 ± 0.004	0.810 ± 0.001	0.808 ± 0.003	0.728 ± 0.002	0.770 ± 0.001	0.819 ± 0.002	0.848 ± 0.001	0.840 ± 0.002
ED_dispo cls	0.734 ± 0.003	0.710 ± 0.006	0.720 ± 0.003	0.676 ± 0.003	0.702 ± 0.005	0.732 ± 0.002	0.750 ± 0.002	0.750 ± 0.002
Potassium reg	0.238 ± 0.017	0.333 ± 0.009	0.326 ± 0.005	0.216 ± 0.010	0.192 ± 0.019	0.308 ± 0.010	0.389 ± 0.005	0.386 ± 0.006
Calcium reg	0.203 ± 0.010	0.181 ± 0.011	0.180 ± 0.009	0.096 ± 0.003	0.130 ± 0.021	0.189 ± 0.007	0.231 ± 0.015	0.265 ± 0.010
eGFR reg	0.573 ± 0.008	0.515 ± 0.008	0.541 ± 0.007	0.448 ± 0.008	0.488 ± 0.017	0.585 ± 0.006	0.596 ± 0.003	0.613 ± 0.003
Glucose reg	0.311 ± 0.014	0.254 ± 0.010	0.279 ± 0.008	0.205 ± 0.010	0.197 ± 0.022	0.295 ± 0.007	0.303 ± 0.010	0.322 ± 0.010
Hemoglobin reg	0.416 ± 0.007	0.335 ± 0.009	0.340 ± 0.007	0.281 ± 0.007	0.291 ± 0.017	0.409 ± 0.014	0.413 ± 0.009	0.423 ± 0.007
Albumin reg	0.428 ± 0.005	0.327 ± 0.005	0.365 ± 0.010	0.310 ± 0.007	0.283 ± 0.011	0.423 ± 0.003	0.427 ± 0.006	0.450 ± 0.006
BUN reg	0.386 ± 0.007	0.327 ± 0.008	0.354 ± 0.013	0.263 ± 0.017	0.272 ± 0.011	0.384 ± 0.005	0.398 ± 0.008	0.410 ± 0.004
Sodium reg	0.240 ± 0.016	0.200 ± 0.017	0.204 ± 0.009	0.168 ± 0.014	0.158 ± 0.010	0.220 ± 0.009	0.265 ± 0.004	0.272 ± 0.006
I48 cls	0.814 ± 0.006	0.827 ± 0.010	0.818 ± 0.008	0.782 ± 0.003	0.836 ± 0.007	0.841 ± 0.008	0.840 ± 0.010	0.857 ± 0.005
I50 cls	0.834 ± 0.005	0.813 ± 0.007	0.815 ± 0.015	0.782 ± 0.003	0.831 ± 0.007	0.832 ± 0.007	0.842 ± 0.002	0.839 ± 0.004
{Z95.0, Z95.2, Z95.3} cls	0.776 ± 0.010	0.796 ± 0.012	0.788 ± 0.006	0.755 ± 0.007	0.817 ± 0.017	0.837 ± 0.007	0.815 ± 0.012	0.852 ± 0.009
N18 cls	0.782 ± 0.015	0.727 ± 0.009	0.742 ± 0.011	0.728 ± 0.007	0.733 ± 0.006	0.764 ± 0.005	0.778 ± 0.005	0.782 ± 0.003
{E10, E11} cls	0.717 ± 0.008	0.666 ± 0.006	0.670 ± 0.009	0.638 ± 0.011	0.667 ± 0.004	0.695 ± 0.008	0.721 ± 0.006	0.720 ± 0.006
{G47.3, E66.2} cls	0.699 ± 0.018	0.676 ± 0.012	0.692 ± 0.012	0.661 ± 0.008	0.669 ± 0.015	0.710 ± 0.014	0.715 ± 0.010	0.715 ± 0.017
{D50, D51, D52, D53, D64} cls	0.677 ± 0.010	0.619 ± 0.007	0.620 ± 0.008	0.637 ± 0.005	0.616 ± 0.007	0.652 ± 0.005	0.668 ± 0.007	0.669 ± 0.011
{J96, R09.02, R06.0} cls	0.691 ± 0.005	0.653 ± 0.007	0.654 ± 0.008	0.647 ± 0.007	0.659 ± 0.004	0.683 ± 0.006	0.684 ± 0.005	0.685 ± 0.002
{I26, I80, I82} cls	0.630 ± 0.015	0.636 ± 0.011	0.634 ± 0.022	0.614 ± 0.014	0.579 ± 0.016	0.611 ± 0.008	0.663 ± 0.014	0.685 ± 0.015

Table S7. **Model comparison of ECG-only *val* performance on all 20 benchmark tasks when using 50% of the train visits.** Tasks are grouped into two demographics tasks, ED disposition classification, eight laboratory value regression tasks, and nine prior ICD-10 diagnosis classification tasks. Performance on the *val* set when using 50% of the train set visits. Training visit are sampled with replacement and both downstream training and hyperparameter selection are repeated five times, results are mean \pm std across repetitions. Regression tasks are evaluated using Pearson correlation (\uparrow), and classification tasks using AUROC (\uparrow). **Bold red** marks the best mean value in each row.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.694 \pm 0.001	0.685 \pm 0.003	0.699 \pm 0.004	0.554 \pm 0.003	0.672 \pm 0.009	0.775 \pm 0.002	0.776 \pm 0.002	0.780 \pm 0.002
Sex cls	0.674 \pm 0.005	0.766 \pm 0.006	0.757 \pm 0.006	0.673 \pm 0.006	0.725 \pm 0.002	0.779 \pm 0.003	0.817 \pm 0.002	0.811 \pm 0.006
ED_dispo cls	0.728 \pm 0.002	0.699 \pm 0.008	0.704 \pm 0.003	0.653 \pm 0.003	0.691 \pm 0.005	0.720 \pm 0.004	0.733 \pm 0.004	0.737 \pm 0.006
Potassium reg	0.268 \pm 0.009	0.302 \pm 0.017	0.279 \pm 0.008	0.179 \pm 0.005	0.178 \pm 0.022	0.294 \pm 0.026	0.370 \pm 0.010	0.377 \pm 0.005
Calcium reg	0.150 \pm 0.016	0.133 \pm 0.016	0.109 \pm 0.012	0.088 \pm 0.012	0.077 \pm 0.025	0.142 \pm 0.011	0.199 \pm 0.016	0.195 \pm 0.005
eGFR reg	0.540 \pm 0.005	0.493 \pm 0.004	0.500 \pm 0.006	0.352 \pm 0.009	0.459 \pm 0.019	0.559 \pm 0.007	0.581 \pm 0.003	0.580 \pm 0.005
Glucose reg	0.297 \pm 0.002	0.185 \pm 0.023	0.188 \pm 0.026	0.152 \pm 0.013	0.160 \pm 0.025	0.233 \pm 0.025	0.270 \pm 0.011	0.249 \pm 0.011
Hemoglobin reg	0.350 \pm 0.006	0.258 \pm 0.015	0.254 \pm 0.012	0.199 \pm 0.008	0.234 \pm 0.018	0.304 \pm 0.061	0.356 \pm 0.013	0.366 \pm 0.014
Albumin reg	0.433 \pm 0.006	0.271 \pm 0.009	0.289 \pm 0.009	0.272 \pm 0.013	0.241 \pm 0.019	0.377 \pm 0.042	0.406 \pm 0.010	0.410 \pm 0.003
BUN reg	0.382 \pm 0.008	0.310 \pm 0.023	0.312 \pm 0.012	0.225 \pm 0.010	0.257 \pm 0.018	0.359 \pm 0.012	0.374 \pm 0.009	0.387 \pm 0.012
Sodium reg	0.286 \pm 0.017	0.176 \pm 0.009	0.156 \pm 0.012	0.127 \pm 0.011	0.148 \pm 0.024	0.211 \pm 0.013	0.241 \pm 0.007	0.255 \pm 0.021
I48 cls	0.825 \pm 0.005	0.813 \pm 0.010	0.795 \pm 0.008	0.752 \pm 0.006	0.840 \pm 0.005	0.854 \pm 0.008	0.842 \pm 0.004	0.864 \pm 0.005
I50 cls	0.845 \pm 0.004	0.810 \pm 0.009	0.796 \pm 0.009	0.773 \pm 0.006	0.823 \pm 0.009	0.819 \pm 0.009	0.837 \pm 0.006	0.837 \pm 0.007
{Z95.0, Z95.2, Z95.3} cls	0.820 \pm 0.007	0.766 \pm 0.019	0.755 \pm 0.019	0.726 \pm 0.009	0.820 \pm 0.006	0.830 \pm 0.006	0.817 \pm 0.010	0.844 \pm 0.008
N18 cls	0.797 \pm 0.003	0.709 \pm 0.021	0.710 \pm 0.010	0.707 \pm 0.010	0.731 \pm 0.007	0.758 \pm 0.011	0.773 \pm 0.008	0.782 \pm 0.004
{E10, E11} cls	0.709 \pm 0.005	0.623 \pm 0.022	0.631 \pm 0.009	0.620 \pm 0.007	0.642 \pm 0.009	0.681 \pm 0.007	0.694 \pm 0.010	0.694 \pm 0.008
{G47.3, E66.2} cls	0.673 \pm 0.013	0.643 \pm 0.008	0.620 \pm 0.020	0.632 \pm 0.014	0.653 \pm 0.011	0.669 \pm 0.005	0.684 \pm 0.009	0.680 \pm 0.013
{D50, D51, D52, D53, D64} cls	0.671 \pm 0.003	0.592 \pm 0.007	0.586 \pm 0.012	0.615 \pm 0.011	0.602 \pm 0.010	0.652 \pm 0.008	0.661 \pm 0.006	0.673 \pm 0.011
{J96, R09.02, R06.0} cls	0.676 \pm 0.002	0.639 \pm 0.010	0.636 \pm 0.007	0.633 \pm 0.007	0.645 \pm 0.005	0.657 \pm 0.015	0.673 \pm 0.007	0.676 \pm 0.007
{I26, I80, I82} cls	0.670 \pm 0.018	0.603 \pm 0.016	0.580 \pm 0.013	0.608 \pm 0.010	0.593 \pm 0.017	0.611 \pm 0.021	0.654 \pm 0.018	0.678 \pm 0.017

Table S8. **Model comparison of PPG-only *val* performance on all 20 benchmark tasks when using 50% of the train visits.** Identical format and content as Table S7, but for *PPG-only*. Performance on the *val* set when using 50% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.718 \pm 0.013	0.701 \pm 0.004	0.716 \pm 0.002	0.674 \pm 0.003	0.604 \pm 0.068	0.655 \pm 0.083	0.768 \pm 0.003	0.793 \pm 0.002
Sex cls	0.740 \pm 0.004	0.729 \pm 0.002	0.726 \pm 0.004	0.692 \pm 0.001	0.711 \pm 0.008	0.734 \pm 0.006	0.764 \pm 0.006	0.765 \pm 0.004
ED_dispo cls	0.704 \pm 0.004	0.687 \pm 0.004	0.682 \pm 0.002	0.670 \pm 0.004	0.664 \pm 0.003	0.677 \pm 0.008	0.703 \pm 0.002	0.721 \pm 0.004
Potassium reg	0.136 \pm 0.008	0.173 \pm 0.010	0.167 \pm 0.018	0.172 \pm 0.007	0.095 \pm 0.036	0.136 \pm 0.019	0.169 \pm 0.014	0.175 \pm 0.008
Calcium reg	0.133 \pm 0.018	0.125 \pm 0.016	0.157 \pm 0.010	0.101 \pm 0.006	0.100 \pm 0.015	0.119 \pm 0.026	0.171 \pm 0.015	0.188 \pm 0.018
eGFR reg	0.467 \pm 0.009	0.425 \pm 0.005	0.458 \pm 0.006	0.424 \pm 0.002	0.322 \pm 0.020	0.381 \pm 0.041	0.503 \pm 0.005	0.552 \pm 0.006
Glucose reg	0.245 \pm 0.042	0.242 \pm 0.020	0.258 \pm 0.029	0.245 \pm 0.008	0.160 \pm 0.025	0.220 \pm 0.034	0.273 \pm 0.012	0.272 \pm 0.007
Hemoglobin reg	0.294 \pm 0.014	0.302 \pm 0.006	0.295 \pm 0.017	0.314 \pm 0.016	0.182 \pm 0.019	0.299 \pm 0.017	0.341 \pm 0.007	0.363 \pm 0.015
Albumin reg	0.302 \pm 0.011	0.269 \pm 0.014	0.309 \pm 0.007	0.304 \pm 0.012	0.173 \pm 0.027	0.264 \pm 0.025	0.326 \pm 0.013	0.384 \pm 0.008
BUN reg	0.267 \pm 0.019	0.239 \pm 0.012	0.282 \pm 0.004	0.281 \pm 0.010	0.162 \pm 0.048	0.234 \pm 0.050	0.311 \pm 0.005	0.354 \pm 0.016
Sodium reg	0.176 \pm 0.027	0.197 \pm 0.021	0.218 \pm 0.018	0.220 \pm 0.011	0.122 \pm 0.039	0.151 \pm 0.022	0.212 \pm 0.008	0.229 \pm 0.006
I48 cls	0.766 \pm 0.016	0.768 \pm 0.010	0.752 \pm 0.006	0.771 \pm 0.008	0.741 \pm 0.012	0.770 \pm 0.007	0.767 \pm 0.010	0.805 \pm 0.012
I50 cls	0.727 \pm 0.005	0.721 \pm 0.011	0.713 \pm 0.010	0.735 \pm 0.006	0.684 \pm 0.009	0.724 \pm 0.003	0.739 \pm 0.007	0.774 \pm 0.007
{Z95.0, Z95.2, Z95.3} cls	0.667 \pm 0.015	0.665 \pm 0.008	0.646 \pm 0.011	0.658 \pm 0.011	0.658 \pm 0.011	0.685 \pm 0.010	0.684 \pm 0.018	0.725 \pm 0.006
N18 cls	0.694 \pm 0.008	0.651 \pm 0.004	0.659 \pm 0.017	0.667 \pm 0.007	0.657 \pm 0.015	0.671 \pm 0.008	0.702 \pm 0.017	0.745 \pm 0.004
{E10, E11} cls	0.679 \pm 0.011	0.655 \pm 0.008	0.659 \pm 0.008	0.631 \pm 0.007	0.643 \pm 0.006	0.669 \pm 0.008	0.685 \pm 0.007	0.700 \pm 0.004
{G47.3, E66.2} cls	0.649 \pm 0.017	0.656 \pm 0.019	0.662 \pm 0.007	0.653 \pm 0.014	0.610 \pm 0.006	0.648 \pm 0.023	0.665 \pm 0.007	0.690 \pm 0.012
{D50, D51, D52, D53, D64} cls	0.627 \pm 0.006	0.605 \pm 0.016	0.594 \pm 0.013	0.609 \pm 0.006	0.607 \pm 0.013	0.582 \pm 0.010	0.624 \pm 0.009	0.647 \pm 0.012
{J96, R09.02, R06.0} cls	0.657 \pm 0.004	0.639 \pm 0.006	0.634 \pm 0.009	0.646 \pm 0.013	0.612 \pm 0.011	0.620 \pm 0.008	0.629 \pm 0.008	0.660 \pm 0.007
{I26, I80, I82} cls	0.621 \pm 0.015	0.636 \pm 0.012	0.595 \pm 0.019	0.626 \pm 0.011	0.566 \pm 0.016	0.593 \pm 0.015	0.612 \pm 0.005	0.646 \pm 0.014

Table S9. **Model comparison of ECG + PPG *val* performance on all 20 benchmark tasks when using 50% of the train visits.** Identical format and content as Table S7, but for *ECG + PPG*. Performance on the *val* set when using 50% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.797 ± 0.003	0.768 ± 0.003	0.774 ± 0.002	0.687 ± 0.005	0.730 ± 0.011	0.818 ± 0.003	0.823 ± 0.002	0.839 ± 0.002
Sex cls	0.767 ± 0.003	0.798 ± 0.004	0.798 ± 0.005	0.727 ± 0.004	0.765 ± 0.002	0.815 ± 0.003	0.842 ± 0.002	0.835 ± 0.002
ED_dispo cls	0.731 ± 0.005	0.708 ± 0.004	0.715 ± 0.005	0.673 ± 0.005	0.695 ± 0.005	0.724 ± 0.003	0.736 ± 0.004	0.745 ± 0.006
Potassium reg	0.223 ± 0.017	0.315 ± 0.012	0.297 ± 0.007	0.199 ± 0.007	0.189 ± 0.012	0.269 ± 0.035	0.384 ± 0.012	0.376 ± 0.006
Calcium reg	0.192 ± 0.020	0.165 ± 0.017	0.165 ± 0.014	0.104 ± 0.015	0.100 ± 0.022	0.164 ± 0.008	0.220 ± 0.017	0.254 ± 0.010
eGFR reg	0.568 ± 0.004	0.508 ± 0.008	0.529 ± 0.004	0.423 ± 0.011	0.467 ± 0.025	0.564 ± 0.008	0.588 ± 0.003	0.598 ± 0.004
Glucose reg	0.298 ± 0.029	0.254 ± 0.020	0.264 ± 0.016	0.202 ± 0.011	0.218 ± 0.035	0.281 ± 0.014	0.313 ± 0.014	0.300 ± 0.009
Hemoglobin reg	0.391 ± 0.015	0.324 ± 0.012	0.323 ± 0.010	0.265 ± 0.012	0.270 ± 0.014	0.386 ± 0.008	0.396 ± 0.012	0.410 ± 0.011
Albumin reg	0.412 ± 0.017	0.307 ± 0.004	0.339 ± 0.006	0.295 ± 0.010	0.272 ± 0.014	0.381 ± 0.034	0.409 ± 0.010	0.429 ± 0.005
BUN reg	0.372 ± 0.021	0.318 ± 0.027	0.335 ± 0.011	0.264 ± 0.011	0.259 ± 0.026	0.354 ± 0.007	0.388 ± 0.013	0.406 ± 0.011
Sodium reg	0.230 ± 0.032	0.195 ± 0.016	0.186 ± 0.017	0.160 ± 0.010	0.161 ± 0.020	0.200 ± 0.019	0.243 ± 0.014	0.269 ± 0.023
I48 cls	0.814 ± 0.014	0.819 ± 0.010	0.806 ± 0.007	0.775 ± 0.007	0.835 ± 0.003	0.848 ± 0.007	0.833 ± 0.005	0.859 ± 0.004
I50 cls	0.827 ± 0.008	0.811 ± 0.010	0.802 ± 0.004	0.787 ± 0.008	0.824 ± 0.010	0.827 ± 0.005	0.837 ± 0.005	0.837 ± 0.008
{Z95.0, Z95.2, Z95.3} cls	0.777 ± 0.027	0.778 ± 0.005	0.762 ± 0.013	0.742 ± 0.012	0.808 ± 0.011	0.822 ± 0.009	0.800 ± 0.008	0.840 ± 0.009
N18 cls	0.785 ± 0.004	0.722 ± 0.016	0.728 ± 0.012	0.723 ± 0.011	0.732 ± 0.012	0.764 ± 0.011	0.772 ± 0.011	0.785 ± 0.009
{E10, E11} cls	0.719 ± 0.018	0.664 ± 0.008	0.674 ± 0.003	0.634 ± 0.006	0.669 ± 0.009	0.700 ± 0.005	0.719 ± 0.007	0.719 ± 0.008
{G47.3, E66.2} cls	0.688 ± 0.028	0.677 ± 0.011	0.683 ± 0.021	0.647 ± 0.016	0.659 ± 0.012	0.693 ± 0.016	0.706 ± 0.008	0.715 ± 0.007
{D50, D51, D52, D53, D64} cls	0.668 ± 0.011	0.604 ± 0.009	0.607 ± 0.013	0.627 ± 0.007	0.603 ± 0.010	0.647 ± 0.008	0.663 ± 0.005	0.665 ± 0.008
{J96, R09.02, R06.0} cls	0.685 ± 0.005	0.657 ± 0.009	0.651 ± 0.010	0.650 ± 0.008	0.645 ± 0.006	0.670 ± 0.013	0.674 ± 0.007	0.685 ± 0.006
{I26, I80, I82} cls	0.646 ± 0.023	0.623 ± 0.015	0.612 ± 0.017	0.615 ± 0.014	0.580 ± 0.021	0.609 ± 0.018	0.658 ± 0.011	0.690 ± 0.020

Table S10. **Model comparison of ECG-only *val* performance on all 20 benchmark tasks when using 25% of the train visits.** Tasks are grouped into two demographics tasks, ED disposition classification, eight laboratory value regression tasks, and nine prior ICD-10 diagnosis classification tasks. Performance on the *val* set when using 25% of the train set visits. Training visit are sampled with replacement and both downstream training and hyperparameter selection are repeated five times, results are mean ± std across repetitions. Regression tasks are evaluated using Pearson correlation (\uparrow), and classification tasks using AUROC (\uparrow). **Bold red** marks the best mean value in each row.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.689 ± 0.004	0.677 ± 0.005	0.681 ± 0.003	0.538 ± 0.007	0.645 ± 0.026	0.762 ± 0.013	0.769 ± 0.002	0.775 ± 0.004
Sex cls	0.665 ± 0.006	0.753 ± 0.003	0.737 ± 0.004	0.668 ± 0.005	0.713 ± 0.004	0.764 ± 0.004	0.806 ± 0.003	0.802 ± 0.004
ED_dispo cls	0.725 ± 0.002	0.691 ± 0.004	0.690 ± 0.007	0.644 ± 0.003	0.693 ± 0.003	0.711 ± 0.006	0.724 ± 0.002	0.730 ± 0.005
Potassium reg	0.239 ± 0.018	0.291 ± 0.012	0.262 ± 0.018	0.175 ± 0.025	0.187 ± 0.020	0.285 ± 0.008	0.345 ± 0.014	0.362 ± 0.009
Calcium reg	0.141 ± 0.016	0.111 ± 0.011	0.101 ± 0.014	0.069 ± 0.010	0.050 ± 0.039	0.133 ± 0.013	0.183 ± 0.008	0.189 ± 0.012
eGFR reg	0.536 ± 0.005	0.468 ± 0.005	0.474 ± 0.010	0.352 ± 0.010	0.467 ± 0.021	0.554 ± 0.014	0.558 ± 0.007	0.567 ± 0.005
Glucose reg	0.278 ± 0.019	0.161 ± 0.023	0.154 ± 0.019	0.128 ± 0.028	0.148 ± 0.036	0.183 ± 0.044	0.235 ± 0.021	0.241 ± 0.021
Hemoglobin reg	0.352 ± 0.006	0.234 ± 0.012	0.226 ± 0.018	0.180 ± 0.014	0.223 ± 0.022	0.301 ± 0.038	0.335 ± 0.009	0.359 ± 0.019
Albumin reg	0.426 ± 0.006	0.263 ± 0.008	0.273 ± 0.012	0.256 ± 0.008	0.252 ± 0.011	0.343 ± 0.109	0.394 ± 0.008	0.414 ± 0.007
BUN reg	0.368 ± 0.010	0.285 ± 0.014	0.274 ± 0.019	0.215 ± 0.010	0.244 ± 0.061	0.331 ± 0.037	0.352 ± 0.013	0.375 ± 0.010
Sodium reg	0.270 ± 0.008	0.162 ± 0.026	0.125 ± 0.009	0.113 ± 0.014	0.146 ± 0.024	0.217 ± 0.023	0.224 ± 0.019	0.250 ± 0.010
I48 cls	0.822 ± 0.004	0.803 ± 0.012	0.786 ± 0.009	0.742 ± 0.003	0.838 ± 0.004	0.843 ± 0.009	0.837 ± 0.004	0.859 ± 0.001
I50 cls	0.838 ± 0.004	0.807 ± 0.007	0.792 ± 0.009	0.753 ± 0.010	0.819 ± 0.008	0.814 ± 0.012	0.823 ± 0.005	0.833 ± 0.006
{Z95.0, Z95.2, Z95.3} cls	0.817 ± 0.008	0.795 ± 0.011	0.764 ± 0.020	0.719 ± 0.009	0.824 ± 0.009	0.833 ± 0.013	0.815 ± 0.005	0.851 ± 0.011
N18 cls	0.791 ± 0.010	0.707 ± 0.016	0.703 ± 0.013	0.690 ± 0.011	0.709 ± 0.008	0.744 ± 0.016	0.746 ± 0.014	0.764 ± 0.010
{E10, E11} cls	0.699 ± 0.010	0.633 ± 0.007	0.618 ± 0.011	0.619 ± 0.005	0.630 ± 0.010	0.669 ± 0.013	0.676 ± 0.009	0.689 ± 0.009
{G47.3, E66.2} cls	0.648 ± 0.022	0.631 ± 0.015	0.609 ± 0.013	0.612 ± 0.009	0.645 ± 0.014	0.649 ± 0.036	0.660 ± 0.025	0.672 ± 0.018
{D50, D51, D52, D53, D64} cls	0.664 ± 0.009	0.605 ± 0.002	0.591 ± 0.009	0.608 ± 0.013	0.601 ± 0.013	0.638 ± 0.016	0.640 ± 0.002	0.672 ± 0.015
{J96, R09.02, R06.0} cls	0.671 ± 0.006	0.627 ± 0.025	0.625 ± 0.012	0.617 ± 0.012	0.637 ± 0.019	0.649 ± 0.012	0.650 ± 0.010	0.659 ± 0.013
{I26, I80, I82} cls	0.639 ± 0.008	0.564 ± 0.019	0.580 ± 0.025	0.591 ± 0.019	0.568 ± 0.025	0.589 ± 0.022	0.610 ± 0.019	0.653 ± 0.028

Table S11. Model comparison of PPG-only *val* performance on all 20 benchmark tasks when using 25% of the train visits. Identical format and content as Table S10, but for PPG-only. Performance on the *val* set when using 25% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.711 ± 0.007	0.690 ± 0.003	0.704 ± 0.003	0.661 ± 0.003	0.580 ± 0.082	0.696 ± 0.015	0.754 ± 0.003	0.783 ± 0.004
Sex cls	0.729 ± 0.005	0.724 ± 0.004	0.717 ± 0.004	0.688 ± 0.003	0.699 ± 0.006	0.723 ± 0.004	0.754 ± 0.004	0.757 ± 0.003
ED_dispo cls	0.689 ± 0.007	0.678 ± 0.005	0.676 ± 0.007	0.651 ± 0.008	0.660 ± 0.007	0.674 ± 0.009	0.696 ± 0.005	0.720 ± 0.002
Potassium reg	0.092 ± 0.036	0.156 ± 0.018	0.144 ± 0.018	0.148 ± 0.008	0.088 ± 0.024	0.099 ± 0.026	0.146 ± 0.022	0.154 ± 0.019
Calcium reg	0.107 ± 0.021	0.106 ± 0.010	0.130 ± 0.017	0.081 ± 0.020	0.085 ± 0.018	0.079 ± 0.033	0.168 ± 0.017	0.188 ± 0.027
eGFR reg	0.412 ± 0.021	0.401 ± 0.008	0.422 ± 0.013	0.419 ± 0.006	0.318 ± 0.027	0.407 ± 0.015	0.476 ± 0.008	0.539 ± 0.002
Glucose reg	0.204 ± 0.056	0.216 ± 0.020	0.245 ± 0.034	0.234 ± 0.011	0.170 ± 0.034	0.224 ± 0.021	0.239 ± 0.016	0.256 ± 0.022
Hemoglobin reg	0.242 ± 0.035	0.280 ± 0.012	0.268 ± 0.019	0.298 ± 0.012	0.189 ± 0.036	0.260 ± 0.026	0.305 ± 0.012	0.337 ± 0.009
Albumin reg	0.227 ± 0.052	0.250 ± 0.014	0.284 ± 0.016	0.290 ± 0.003	0.150 ± 0.036	0.208 ± 0.036	0.292 ± 0.024	0.357 ± 0.012
BUN reg	0.193 ± 0.027	0.212 ± 0.013	0.233 ± 0.011	0.242 ± 0.010	0.120 ± 0.045	0.196 ± 0.018	0.267 ± 0.019	0.323 ± 0.013
Sodium reg	0.128 ± 0.023	0.175 ± 0.021	0.185 ± 0.007	0.205 ± 0.013	0.091 ± 0.036	0.113 ± 0.043	0.191 ± 0.011	0.221 ± 0.019
I48 cls	0.744 ± 0.010	0.758 ± 0.002	0.754 ± 0.009	0.739 ± 0.012	0.734 ± 0.009	0.771 ± 0.007	0.759 ± 0.009	0.793 ± 0.008
I50 cls	0.703 ± 0.004	0.715 ± 0.009	0.700 ± 0.008	0.710 ± 0.006	0.669 ± 0.012	0.701 ± 0.010	0.714 ± 0.009	0.757 ± 0.012
{Z95.0, Z95.2, Z95.3} cls	0.652 ± 0.020	0.655 ± 0.025	0.629 ± 0.024	0.637 ± 0.012	0.647 ± 0.017	0.672 ± 0.019	0.658 ± 0.013	0.705 ± 0.019
N18 cls	0.674 ± 0.011	0.641 ± 0.014	0.662 ± 0.006	0.656 ± 0.013	0.633 ± 0.014	0.652 ± 0.011	0.686 ± 0.014	0.722 ± 0.009
{E10, E11} cls	0.673 ± 0.012	0.638 ± 0.010	0.639 ± 0.021	0.603 ± 0.010	0.631 ± 0.010	0.667 ± 0.014	0.668 ± 0.009	0.677 ± 0.011
{G47.3, E66.2} cls	0.652 ± 0.022	0.657 ± 0.010	0.640 ± 0.021	0.655 ± 0.005	0.608 ± 0.013	0.655 ± 0.009	0.670 ± 0.024	0.678 ± 0.022
{D50, D51, D52, D53, D64} cls	0.620 ± 0.004	0.610 ± 0.014	0.592 ± 0.014	0.603 ± 0.009	0.600 ± 0.021	0.599 ± 0.020	0.626 ± 0.007	0.649 ± 0.015
{J96, R09.02, R06.0} cls	0.636 ± 0.014	0.626 ± 0.003	0.613 ± 0.008	0.633 ± 0.005	0.598 ± 0.010	0.617 ± 0.011	0.621 ± 0.007	0.650 ± 0.008
{I26, I80, I82} cls	0.600 ± 0.016	0.625 ± 0.019	0.598 ± 0.022	0.618 ± 0.006	0.558 ± 0.014	0.602 ± 0.011	0.618 ± 0.011	0.646 ± 0.020

Table S12. Model comparison of ECG + PPG *val* performance on all 20 benchmark tasks when using 25% of the train visits. Identical format and content as Table S10, but for ECG + PPG. Performance on the *val* set when using 25% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.789 ± 0.003	0.759 ± 0.004	0.762 ± 0.003	0.669 ± 0.002	0.696 ± 0.026	0.798 ± 0.028	0.816 ± 0.002	0.833 ± 0.002
Sex cls	0.757 ± 0.006	0.785 ± 0.002	0.781 ± 0.002	0.720 ± 0.004	0.748 ± 0.004	0.800 ± 0.003	0.829 ± 0.002	0.826 ± 0.004
ED_dispo cls	0.722 ± 0.007	0.697 ± 0.005	0.701 ± 0.004	0.660 ± 0.003	0.697 ± 0.003	0.719 ± 0.005	0.728 ± 0.003	0.737 ± 0.004
Potassium reg	0.157 ± 0.034	0.297 ± 0.007	0.274 ± 0.014	0.187 ± 0.030	0.191 ± 0.029	0.258 ± 0.014	0.360 ± 0.014	0.361 ± 0.004
Calcium reg	0.167 ± 0.012	0.141 ± 0.011	0.144 ± 0.023	0.082 ± 0.012	0.074 ± 0.044	0.152 ± 0.012	0.203 ± 0.012	0.248 ± 0.012
eGFR reg	0.523 ± 0.014	0.492 ± 0.005	0.508 ± 0.013	0.418 ± 0.007	0.469 ± 0.031	0.549 ± 0.025	0.570 ± 0.004	0.586 ± 0.002
Glucose reg	0.243 ± 0.045	0.238 ± 0.018	0.232 ± 0.016	0.180 ± 0.031	0.186 ± 0.040	0.222 ± 0.039	0.283 ± 0.023	0.288 ± 0.016
Hemoglobin reg	0.351 ± 0.031	0.298 ± 0.013	0.289 ± 0.020	0.259 ± 0.019	0.238 ± 0.028	0.334 ± 0.046	0.376 ± 0.014	0.396 ± 0.013
Albumin reg	0.350 ± 0.036	0.294 ± 0.004	0.317 ± 0.012	0.287 ± 0.003	0.254 ± 0.018	0.374 ± 0.042	0.398 ± 0.006	0.425 ± 0.010
BUN reg	0.303 ± 0.026	0.288 ± 0.015	0.289 ± 0.013	0.254 ± 0.009	0.238 ± 0.070	0.293 ± 0.057	0.359 ± 0.016	0.385 ± 0.012
Sodium reg	0.180 ± 0.029	0.177 ± 0.023	0.155 ± 0.012	0.147 ± 0.014	0.154 ± 0.031	0.202 ± 0.025	0.230 ± 0.017	0.257 ± 0.011
I48 cls	0.797 ± 0.012	0.808 ± 0.013	0.791 ± 0.010	0.768 ± 0.003	0.834 ± 0.005	0.837 ± 0.007	0.832 ± 0.006	0.854 ± 0.003
I50 cls	0.813 ± 0.013	0.807 ± 0.008	0.799 ± 0.013	0.763 ± 0.006	0.819 ± 0.009	0.817 ± 0.009	0.826 ± 0.007	0.828 ± 0.007
{Z95.0, Z95.2, Z95.3} cls	0.766 ± 0.025	0.791 ± 0.024	0.752 ± 0.015	0.731 ± 0.011	0.815 ± 0.017	0.822 ± 0.012	0.813 ± 0.013	0.845 ± 0.012
N18 cls	0.766 ± 0.019	0.710 ± 0.023	0.719 ± 0.011	0.704 ± 0.008	0.709 ± 0.009	0.741 ± 0.009	0.748 ± 0.014	0.762 ± 0.010
{E10, E11} cls	0.715 ± 0.015	0.658 ± 0.012	0.655 ± 0.008	0.635 ± 0.005	0.648 ± 0.011	0.688 ± 0.009	0.701 ± 0.007	0.703 ± 0.010
{G47.3, E66.2} cls	0.689 ± 0.022	0.656 ± 0.014	0.646 ± 0.021	0.636 ± 0.013	0.658 ± 0.021	0.676 ± 0.028	0.685 ± 0.027	0.700 ± 0.008
{D50, D51, D52, D53, D64} cls	0.665 ± 0.006	0.619 ± 0.001	0.606 ± 0.010	0.622 ± 0.016	0.604 ± 0.012	0.638 ± 0.010	0.647 ± 0.010	0.663 ± 0.011
{J96, R09.02, R06.0} cls	0.667 ± 0.018	0.637 ± 0.023	0.631 ± 0.009	0.632 ± 0.012	0.636 ± 0.025	0.655 ± 0.010	0.658 ± 0.013	0.662 ± 0.012
{I26, I80, I82} cls	0.612 ± 0.019	0.589 ± 0.027	0.602 ± 0.020	0.602 ± 0.016	0.555 ± 0.027	0.590 ± 0.021	0.617 ± 0.014	0.657 ± 0.022

Table S13. **Model comparison of ECG-only *val* performance on all 20 benchmark tasks when using 10% of the train visits.** Tasks are grouped into two demographics tasks, ED disposition classification, eight laboratory value regression tasks, and nine prior ICD-10 diagnosis classification tasks. Performance on the *val* set when using 10% of the train set visits. Training visit are sampled with replacement and both downstream training and hyperparameter selection are repeated five times, results are mean \pm std across repetitions. Regression tasks are evaluated using Pearson correlation (\uparrow), and classification tasks using AUROC (\uparrow). **Bold red** marks the best mean value in each row.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.681 \pm 0.003	0.663 \pm 0.002	0.657 \pm 0.004	0.517 \pm 0.008	0.521 \pm 0.139	0.726 \pm 0.042	0.755 \pm 0.004	0.761 \pm 0.003
Sex cls	0.656 \pm 0.011	0.735 \pm 0.006	0.715 \pm 0.004	0.665 \pm 0.009	0.706 \pm 0.004	0.745 \pm 0.007	0.791 \pm 0.004	0.788 \pm 0.003
ED_dispo cls	0.712 \pm 0.011	0.687 \pm 0.005	0.677 \pm 0.014	0.638 \pm 0.007	0.690 \pm 0.010	0.708 \pm 0.015	0.718 \pm 0.004	0.726 \pm 0.009
Potassium reg	0.220 \pm 0.022	0.249 \pm 0.012	0.229 \pm 0.014	0.141 \pm 0.020	0.124 \pm 0.064	0.232 \pm 0.021	0.310 \pm 0.014	0.329 \pm 0.011
Calcium reg	0.125 \pm 0.013	0.094 \pm 0.012	0.100 \pm 0.019	0.062 \pm 0.016	0.060 \pm 0.010	0.111 \pm 0.005	0.163 \pm 0.006	0.163 \pm 0.014
eGFR reg	0.519 \pm 0.010	0.440 \pm 0.015	0.440 \pm 0.007	0.333 \pm 0.016	0.423 \pm 0.038	0.483 \pm 0.060	0.541 \pm 0.005	0.558 \pm 0.011
Glucose reg	0.235 \pm 0.070	0.150 \pm 0.028	0.131 \pm 0.027	0.117 \pm 0.014	0.144 \pm 0.057	0.195 \pm 0.038	0.211 \pm 0.019	0.229 \pm 0.025
Hemoglobin reg	0.313 \pm 0.018	0.196 \pm 0.007	0.225 \pm 0.019	0.154 \pm 0.019	0.169 \pm 0.061	0.270 \pm 0.044	0.310 \pm 0.009	0.350 \pm 0.010
Albumin reg	0.405 \pm 0.011	0.238 \pm 0.033	0.242 \pm 0.021	0.218 \pm 0.016	0.206 \pm 0.046	0.315 \pm 0.048	0.350 \pm 0.016	0.377 \pm 0.012
BUN reg	0.351 \pm 0.014	0.274 \pm 0.017	0.282 \pm 0.022	0.193 \pm 0.026	0.246 \pm 0.060	0.322 \pm 0.044	0.325 \pm 0.024	0.342 \pm 0.040
Sodium reg	0.251 \pm 0.024	0.132 \pm 0.011	0.125 \pm 0.015	0.111 \pm 0.011	0.134 \pm 0.063	0.148 \pm 0.055	0.195 \pm 0.023	0.233 \pm 0.007
I48 cls	0.797 \pm 0.022	0.785 \pm 0.007	0.774 \pm 0.010	0.721 \pm 0.017	0.833 \pm 0.008	0.828 \pm 0.010	0.823 \pm 0.005	0.847 \pm 0.005
I50 cls	0.818 \pm 0.020	0.793 \pm 0.008	0.766 \pm 0.016	0.739 \pm 0.009	0.817 \pm 0.008	0.809 \pm 0.011	0.814 \pm 0.009	0.825 \pm 0.005
{Z95.0, Z95.2, Z95.3} cls	0.789 \pm 0.028	0.748 \pm 0.038	0.733 \pm 0.028	0.700 \pm 0.021	0.801 \pm 0.013	0.821 \pm 0.018	0.787 \pm 0.033	0.838 \pm 0.025
N18 cls	0.754 \pm 0.026	0.689 \pm 0.012	0.683 \pm 0.024	0.667 \pm 0.012	0.692 \pm 0.014	0.725 \pm 0.014	0.731 \pm 0.011	0.750 \pm 0.010
{E10, E11} cls	0.683 \pm 0.017	0.612 \pm 0.015	0.592 \pm 0.011	0.604 \pm 0.013	0.614 \pm 0.010	0.647 \pm 0.005	0.656 \pm 0.009	0.675 \pm 0.009
{G47.3, E66.2} cls	0.621 \pm 0.011	0.613 \pm 0.028	0.595 \pm 0.024	0.585 \pm 0.006	0.630 \pm 0.021	0.621 \pm 0.024	0.628 \pm 0.015	0.662 \pm 0.012
{D50, D51, D52, D53, D64} cls	0.652 \pm 0.019	0.590 \pm 0.029	0.594 \pm 0.008	0.571 \pm 0.009	0.591 \pm 0.004	0.634 \pm 0.008	0.642 \pm 0.009	0.669 \pm 0.013
{J96, R09.02, R06.0} cls	0.667 \pm 0.002	0.603 \pm 0.017	0.611 \pm 0.017	0.600 \pm 0.015	0.631 \pm 0.009	0.638 \pm 0.010	0.640 \pm 0.007	0.645 \pm 0.014
{I26, I80, I82} cls	0.610 \pm 0.021	0.557 \pm 0.035	0.562 \pm 0.020	0.568 \pm 0.022	0.553 \pm 0.032	0.568 \pm 0.041	0.573 \pm 0.016	0.624 \pm 0.018

Table S14. **Model comparison of PPG-only *val* performance on all 20 benchmark tasks when using 10% of the train visits.** Identical format and content as Table S13, but for *PPG-only*. Performance on the *val* set when using 10% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.633 \pm 0.035	0.668 \pm 0.006	0.673 \pm 0.006	0.649 \pm 0.003	0.549 \pm 0.049	0.636 \pm 0.023	0.726 \pm 0.003	0.763 \pm 0.005
Sex cls	0.716 \pm 0.004	0.715 \pm 0.005	0.698 \pm 0.004	0.682 \pm 0.006	0.673 \pm 0.013	0.704 \pm 0.008	0.734 \pm 0.004	0.744 \pm 0.006
ED_dispo cls	0.676 \pm 0.007	0.667 \pm 0.006	0.657 \pm 0.005	0.632 \pm 0.007	0.635 \pm 0.008	0.658 \pm 0.005	0.672 \pm 0.007	0.712 \pm 0.005
Potassium reg	0.064 \pm 0.023	0.124 \pm 0.026	0.122 \pm 0.026	0.129 \pm 0.020	0.048 \pm 0.028	0.103 \pm 0.025	0.112 \pm 0.020	0.127 \pm 0.017
Calcium reg	0.076 \pm 0.020	0.085 \pm 0.027	0.076 \pm 0.028	0.072 \pm 0.006	0.030 \pm 0.010	0.075 \pm 0.014	0.118 \pm 0.009	0.160 \pm 0.013
eGFR reg	0.207 \pm 0.029	0.374 \pm 0.008	0.395 \pm 0.019	0.401 \pm 0.004	0.214 \pm 0.097	0.318 \pm 0.082	0.425 \pm 0.018	0.510 \pm 0.015
Glucose reg	0.095 \pm 0.061	0.205 \pm 0.023	0.202 \pm 0.016	0.222 \pm 0.019	0.063 \pm 0.060	0.226 \pm 0.026	0.203 \pm 0.026	0.234 \pm 0.015
Hemoglobin reg	0.173 \pm 0.055	0.260 \pm 0.018	0.244 \pm 0.024	0.274 \pm 0.014	0.142 \pm 0.044	0.261 \pm 0.023	0.259 \pm 0.014	0.303 \pm 0.012
Albumin reg	0.200 \pm 0.039	0.233 \pm 0.013	0.259 \pm 0.017	0.278 \pm 0.013	0.129 \pm 0.033	0.216 \pm 0.035	0.269 \pm 0.025	0.338 \pm 0.019
BUN reg	0.150 \pm 0.038	0.198 \pm 0.022	0.215 \pm 0.027	0.227 \pm 0.019	0.139 \pm 0.020	0.215 \pm 0.015	0.237 \pm 0.026	0.304 \pm 0.021
Sodium reg	0.076 \pm 0.031	0.168 \pm 0.023	0.161 \pm 0.013	0.177 \pm 0.013	0.036 \pm 0.049	0.118 \pm 0.035	0.154 \pm 0.011	0.196 \pm 0.014
I48 cls	0.735 \pm 0.022	0.736 \pm 0.008	0.752 \pm 0.003	0.711 \pm 0.020	0.721 \pm 0.007	0.758 \pm 0.008	0.755 \pm 0.008	0.789 \pm 0.009
I50 cls	0.695 \pm 0.010	0.687 \pm 0.017	0.677 \pm 0.017	0.679 \pm 0.017	0.648 \pm 0.027	0.679 \pm 0.021	0.688 \pm 0.013	0.722 \pm 0.013
{Z95.0, Z95.2, Z95.3} cls	0.634 \pm 0.011	0.634 \pm 0.021	0.596 \pm 0.019	0.609 \pm 0.013	0.615 \pm 0.012	0.657 \pm 0.017	0.638 \pm 0.017	0.683 \pm 0.022
N18 cls	0.653 \pm 0.020	0.634 \pm 0.010	0.634 \pm 0.013	0.627 \pm 0.020	0.605 \pm 0.009	0.635 \pm 0.011	0.655 \pm 0.010	0.707 \pm 0.012
{E10, E11} cls	0.658 \pm 0.011	0.634 \pm 0.015	0.624 \pm 0.012	0.584 \pm 0.004	0.607 \pm 0.014	0.638 \pm 0.008	0.650 \pm 0.012	0.662 \pm 0.011
{G47.3, E66.2} cls	0.642 \pm 0.016	0.622 \pm 0.022	0.605 \pm 0.035	0.625 \pm 0.020	0.602 \pm 0.012	0.621 \pm 0.019	0.640 \pm 0.029	0.662 \pm 0.021
{D50, D51, D52, D53, D64} cls	0.613 \pm 0.017	0.605 \pm 0.021	0.591 \pm 0.022	0.600 \pm 0.007	0.596 \pm 0.013	0.584 \pm 0.011	0.613 \pm 0.015	0.638 \pm 0.006
{J96, R09.02, R06.0} cls	0.631 \pm 0.009	0.615 \pm 0.005	0.610 \pm 0.005	0.622 \pm 0.007	0.586 \pm 0.011	0.605 \pm 0.007	0.608 \pm 0.014	0.630 \pm 0.006
{I26, I80, I82} cls	0.620 \pm 0.021	0.588 \pm 0.031	0.561 \pm 0.050	0.597 \pm 0.016	0.518 \pm 0.011	0.550 \pm 0.042	0.589 \pm 0.031	0.596 \pm 0.040

Table S15. **Model comparison of ECG + PPG *val* performance on all 20 benchmark tasks when using 10% of the train visits.** Identical format and content as Table S13, but for *ECG + PPG*. Performance on the *val* set when using 10% of the train set visits.

Task	Domain Features	MOMENT-base	Chronos-Bolt-small	PaPaGei	D-BETA	ECGFounder	xECG-10min	CSFM-base
Age reg	0.746 ± 0.022	0.747 ± 0.003	0.742 ± 0.005	0.650 ± 0.010	0.528 ± 0.198	0.793 ± 0.004	0.801 ± 0.002	0.823 ± 0.003
Sex cls	0.746 ± 0.003	0.772 ± 0.006	0.755 ± 0.005	0.714 ± 0.009	0.736 ± 0.007	0.782 ± 0.004	0.813 ± 0.004	0.814 ± 0.004
ED_dispo cls	0.712 ± 0.004	0.696 ± 0.007	0.691 ± 0.009	0.652 ± 0.006	0.686 ± 0.015	0.711 ± 0.013	0.718 ± 0.001	0.728 ± 0.006
Potassium reg	0.124 ± 0.053	0.260 ± 0.026	0.226 ± 0.019	0.166 ± 0.021	0.126 ± 0.067	0.199 ± 0.037	0.316 ± 0.027	0.323 ± 0.017
Calcium reg	0.113 ± 0.030	0.115 ± 0.012	0.121 ± 0.010	0.075 ± 0.011	0.087 ± 0.009	0.124 ± 0.014	0.180 ± 0.011	0.210 ± 0.007
eGFR reg	0.334 ± 0.072	0.469 ± 0.011	0.479 ± 0.008	0.381 ± 0.019	0.422 ± 0.031	0.514 ± 0.037	0.547 ± 0.006	0.572 ± 0.010
Glucose reg	0.139 ± 0.068	0.218 ± 0.031	0.193 ± 0.018	0.167 ± 0.008	0.166 ± 0.057	0.202 ± 0.041	0.260 ± 0.015	0.268 ± 0.024
Hemoglobin reg	0.261 ± 0.052	0.254 ± 0.009	0.268 ± 0.009	0.219 ± 0.019	0.170 ± 0.060	0.284 ± 0.081	0.346 ± 0.007	0.383 ± 0.005
Albumin reg	0.302 ± 0.055	0.272 ± 0.023	0.278 ± 0.027	0.250 ± 0.015	0.207 ± 0.048	0.349 ± 0.031	0.354 ± 0.024	0.393 ± 0.011
BUN reg	0.240 ± 0.040	0.283 ± 0.016	0.283 ± 0.030	0.216 ± 0.031	0.244 ± 0.064	0.290 ± 0.043	0.335 ± 0.019	0.352 ± 0.034
Sodium reg	0.111 ± 0.049	0.149 ± 0.002	0.157 ± 0.017	0.135 ± 0.021	0.138 ± 0.061	0.096 ± 0.043	0.195 ± 0.025	0.234 ± 0.015
I48 cls	0.792 ± 0.023	0.796 ± 0.006	0.783 ± 0.011	0.749 ± 0.012	0.824 ± 0.013	0.826 ± 0.010	0.817 ± 0.005	0.837 ± 0.011
I50 cls	0.806 ± 0.012	0.796 ± 0.004	0.769 ± 0.021	0.746 ± 0.011	0.809 ± 0.013	0.808 ± 0.012	0.809 ± 0.016	0.812 ± 0.019
{Z95.0, Z95.2, Z95.3} cls	0.758 ± 0.019	0.756 ± 0.037	0.723 ± 0.028	0.709 ± 0.021	0.785 ± 0.022	0.810 ± 0.018	0.772 ± 0.049	0.828 ± 0.029
N18 cls	0.740 ± 0.020	0.699 ± 0.011	0.689 ± 0.019	0.672 ± 0.013	0.682 ± 0.026	0.722 ± 0.017	0.722 ± 0.015	0.748 ± 0.014
{E10, E11} cls	0.696 ± 0.007	0.642 ± 0.015	0.623 ± 0.011	0.612 ± 0.012	0.627 ± 0.014	0.667 ± 0.007	0.680 ± 0.011	0.688 ± 0.005
{G47.3, E66.2} cls	0.673 ± 0.008	0.643 ± 0.022	0.625 ± 0.022	0.606 ± 0.012	0.640 ± 0.019	0.641 ± 0.025	0.652 ± 0.012	0.687 ± 0.016
{D50, D51, D52, D53, D64} cls	0.651 ± 0.015	0.617 ± 0.020	0.608 ± 0.012	0.584 ± 0.010	0.594 ± 0.013	0.628 ± 0.011	0.640 ± 0.010	0.657 ± 0.011
{J96, R09.02, R06.0} cls	0.667 ± 0.005	0.625 ± 0.017	0.625 ± 0.010	0.613 ± 0.016	0.634 ± 0.012	0.644 ± 0.016	0.643 ± 0.014	0.637 ± 0.020
{I26, I80, I82} cls	0.621 ± 0.015	0.577 ± 0.032	0.559 ± 0.035	0.577 ± 0.028	0.543 ± 0.037	0.563 ± 0.044	0.579 ± 0.023	0.619 ± 0.028

Table S16. **Number of *train*, *val*, and *test* visits for each of the eight laboratory regression tasks.** To ensure unambiguous supervision, only visits with exactly one recorded measurement of the corresponding laboratory test are included.

Task	Train	Val	Test
Potassium reg (“POTASSIUM”)	13,481	1,700	1,683
Calcium reg (“CALCIUM”)	13,464	1,700	1,677
eGFR reg (“EGFR REFIT WITHOUT RACE (2021)”)	13,407	1,694	1,675
Glucose reg (“GLUCOSE”)	13,100	1,658	1,644
Hemoglobin reg (“HEMOGLOBIN (HGB)”)	13,854	1,772	1,724
Albumin reg (“ALBUMIN”)	13,538	1,729	1,707
BUN reg (“BLOOD UREA NITROGEN (BUN)”)	13,096	1,658	1,644
Sodium reg (“SODIUM”)	13,092	1,657	1,643

Table S17. **Positive class prevalence for each of the eleven binary classification tasks.** Percentage of visits labeled as the positive class in the *train*, *val*, and *test* splits. Values are computed at the visit level and are provided to contextualize reported AUROC performance.

Task	Label / Description	Train	Val	Test
Sex cls	“M” vs “F”	52.2%	52.4%	52.5%
ED_dispo cls	“Discharge” vs {“Inpatient”, “Observation”, “ICU”}	50.5%	51.0%	51.2%
I48 cls	Atrial Fibrillation Detection	10.5%	8.5%	7.6%
I50 cls	Heart Failure Detection	9.8%	7.2%	6.5%
{Z95.0, Z95.2, Z95.3} cls	Cardiac Device Detection	3.6%	2.6%	2.4%
N18 cls	Chronic Kidney Disease Detection	10.5%	7.0%	7.4%
{E10, E11} cls	Diabetes Mellitus Detection	15.4%	12.1%	11.7%
{G47.3, E66.2} cls	Sleep-Related Breathing Disorder Detection	9.3%	6.3%	7.3%
{D50, D51, D52, D53, D64} cls	Anemia Detection	16.1%	11.4%	11.1%
{J96, R09.02, R06.0} cls	Respiratory Compromise Detection	23.3%	15.2%	13.1%
{I26, I80, I82} cls	Venous Thromboembolism Detection	8.6%	5.4%	5.5%