

A Signal–Language Foundation Model for Broad-Spectrum Cardiovascular Assessment from Routine Electrocardiography

**Ziqing Yu^{†1}, Yuhui Tao^{†2,3}, Jiayu Huo^{†4}, Lei Pan¹, Zilong Xiao¹, Juecheng Chen¹, Xiao Li¹,
Jianxuan Li¹, You Zhou⁵, Zhixing Li¹, Cong Wang¹, Beijian Zhang¹, Chen Chen¹¹, Hongyang Lu⁶,
Konstantinos Patlatzoglou⁴, Daniel B. Kramer⁷, Jonathan W. Waks⁸, Yangang Su¹, Fu Siong
Ng^{4,9,10}, Shuo Wang^{2,3}, Yixiu Liang^{1,4}, Junbo Ge¹**

1 Department of Cardiology, Zhongshan Hospital of Fudan University, Shanghai Institute of Cardiovascular Diseases, National Clinical Research Centre for Interventional Medicine, Shanghai, China

2 Digital Medical Research Center, School of Basic Medical Sciences, Fudan University, Shanghai, 200032, China.

3 Shanghai Key Laboratory of Medical Imaging Computing and Computer Assisted Intervention, Shanghai, 200032, China.

4 National Heart and Lung Institute, Imperial College London, Hammersmith Hospital, Du Cane Road, London W120HS, UK

5 Department of Cardiology, Shanghai Geriatric Medical Center, Shanghai, China

6 Cardiac Rhythm Management, Medtronic Technology Center, Medtronic (Shanghai) Ltd., Shanghai, China

7 Richard A. and Susan F. Smith Center for Outcomes Research in Cardiology, Beth Israel Deaconess Medical Center, Harvard Medical School, Boston MA USA

8 Harvard-Thorndike Electrophysiology Institute, Beth Israel Deaconess Medical Center, Harvard Medical School, Boston, MA, USA.

9 Department of Cardiology, Imperial College Healthcare NHS Trust, London, United Kingdom

10 Department of Cardiology, Chelsea and Westminster NHS Foundation Trust, London, United

11 Department of Computer Science and Technology, University of Cambridge, Cambridge, UK

†: Equal contribution as first author

Correspondence:

Fu Siong Ng, Email: f.ng@imperial.ac.uk, Shuo Wang, Email: shuowang26@gmail.com, and Yixiu

Liang, Email: liang.yixiu@zs-hospital.sh.cn

Abstract

Background: Electrocardiography (ECG) is a cornerstone of cardiovascular care, yet conventional AI models are often limited to common arrhythmias, lack generalizability across diverse populations, and fail to address rare or diagnostically elusive conditions. We developed ECG Contrastive Language-Image Pre-training (ECGCLIP), a signal-language contrastive learning framework, to address these gaps by leveraging expert-curated ECG-text pairs for pre-training.

Methods: ECGCLIP was pre-trained on 2,837,962 ECG studies from 1,324,856 patients, using self-supervised contrastive learning to align ECG waveforms with expert-generated diagnostic reports. The model was comprehensively evaluated on a held-out internal test set and nine independent external cohorts (comprising around 1.5 million ECGs) across 89 downstream tasks: 45 ECG tasks (arrhythmias, ischemia, conduction, structure, pacing, electrolyte abnormalities), 39 echocardiographic tasks (valvular disease, ventricular function, etc.), and 5 rare cardiac diseases. Performance was quantified using PRAUC, with comparisons to random initialization (Random Init-R18) and Merl-R18 baselines.

Results: ECGCLIP delivered consistent and enhanced diagnostic performance across both internal and external evaluations. On the internal test set, for common clinical conditions such as atrial fibrillation (PRAUC 0.900) and ST-segment elevation myocardial infarction (PRAUC 0.383), the deeper ECGCLIP-R34 architecture significantly outperformed baselines. This robust discriminative ability generalized robustly across all nine independent external cohorts. Notably, substantial improvements were observed in low-prevalence diseases that are not reliably diagnosed via standard ECG. For Ebstein anomaly, constrictive pericarditis, dextrocardia, and cardiac amyloidosis, ECGCLIP-R34 achieved PRAUC values of 0.253, 0.175, 0.121, and 0.201 respectively on the internal test set, representing multi-fold increases over baselines that performed near chance. Additionally, ECGCLIP demonstrated high data efficiency, matching or exceeding full-dataset baseline performance using only 10% of training data. t-SNE visualization and saliency heatmaps confirmed that ECGCLIP learns clinically meaningful, well-separated feature representations, with attention aligned to established electrocardiographic criteria.

Conclusions: ECGCLIP represents a significant advance in cardiac AI, delivering robust diagnostic performance across routine tasks while enabling the opportunistic screening of echocardiographic and

rare conditions. By combining large-scale, expert-curated data with multimodal contrastive learning, it redefines the ECG as a highly scalable clinical gatekeeper, paving the way for more equitable and accessible advanced diagnostics globally.

Keywords: Artificial Intelligence, Electrocardiography, Echocardiography, Rare heart disease

Introduction

Electrocardiography (ECG) remains a fundamental cardiovascular diagnostic tool. Due to its non-invasive nature, cost-effectiveness, wide diagnostic scope, and real-time monitoring capability, it is widely used for the detection and management of various cardiovascular diseases, including arrhythmias, myocardial ischemia, myocardial infarction, and other cardiac conditions¹. Furthermore, the standard ECG harbors latent structural and pathophysiological signals. These include markers for cardiomyopathy, amyloidosis, or cardiac dysfunction. While invisible to conventional visual inspection, these subtle patterns are now detectable by advanced artificial intelligence (AI)². Consequently, foundation models are redefining the ECG from a basic rhythm-monitoring tool into a potential surrogate for structural cardiac assessment.

Despite the rapid integration of AI into cardiovascular medicine^{3,4}, existing AI-ECG approaches remain fundamentally limited. Early supervised models are confined to narrow diagnostic tasks, trained on small, homogeneous cohorts, and often fail to generalize across populations or detect rare or structural diseases^{5,6}. While more recent foundation models have scaled in data volume, they often focus on a restricted set of common arrhythmias or rely on proxy labels, lacking direct alignment with expert clinical interpretations^{7,8}. Crucially, robust performance on echocardiography-derived phenotypes, particularly for rare cardiomyopathies, has rarely been achieved. This gap severely limits their clinical scope and real-world applicability.

The true potential of AI in electrocardiography lies not in replicating known patterns, but in decoding latent structural and functional signals². To unlock this capability, we propose a methodological shift from rigid label-based training to semantic natural language supervision, adapted from the Contrastive Language-Image Pretraining (CLIP) framework. In this study, we introduce ECGCLIP, a signal-language foundation model trained on paired ECG recordings and their corresponding expert-generated diagnostic reports. ECGCLIP leverages the high-fidelity clinical reasoning embedded in human expert narratives to align raw waveforms with complex pathological concepts. We rigorously validated the model across one held-out internal test set and nine independent external cohorts to evaluate its performance across a broad diagnostic spectrum. Specifically, we assess

its capability to handle 45 standard ECG tasks, bridge the modality gap by screening for 39 echocardiographic structural phenotypes, and identify 5 rare cardiomyopathies that require comprehensive clinical adjudication. Our findings position ECGCLIP as a scalable, data-efficient risk stratification tool, aiming to democratize advanced cardiovascular diagnostics and optimize the allocation of advanced imaging resources worldwide.

Methods

Data Resources and Experimental Design

Data Curation and Partitioning Strategy. To ensure robust foundation model learning and rigorous external validation, we curated a multi-institutional dataset comprising over 4.8 million ECG recordings from 10 distinct cohorts (see [Figure 1A](#) and [Table S1](#) for details). The primary dataset for model development originates from Zhongshan Hospital (Zhongshan), encompassing 3,559,279 ECG signal-report pairs collected between 2013 and 2023. To ensure rigorous evaluation and prevent data leakage, we implemented a strict patient-level partitioning strategy. The entire cohort was first randomly divided into a pre-training set and a downstream task set based on an 85:15 patient ratio (see [Table S2](#) for details). Subsequently, internal partitioning was performed within each subset: the pre-training cohort was split into training and validation sets (98:2), while the downstream cohort was further stratified into training, validation, and test sets (70:15:15). Crucially, strict patient independence was maintained across all splits; no individual patient appeared in more than one subset, ensuring that the evaluation reflects generalization capability to unseen subjects. Stratified sampling was applied to preserve the prevalence of rare cardiac conditions, ensuring robust evaluation for long-tail pathologies. For external validation, we incorporated nine independent cohorts: three hospital-based datasets (Shanghai Tenth People's Hospital [Shiyuan], Xiamen Branch of Zhongshan Hospital [Xiamen], and Beth Israel Deaconess Medical Center [BIDMC]) and six public benchmarks (MIMIC-IV-ECG, UK Biobank [UKB], Chapman, PTB-XL, Georgia, and CPSC2018), collectively contributing over 900,000 recordings. All external cohorts were maintained as held-out testing sets to assess model robustness under distribution shifts.

Task Definitions and Label Standardization. We structured the evaluation across three tiers of diagnostic complexity (refer to [Table S3](#) for full task definitions, abbreviations, and hierarchical clustering in the internal development cohort). For all external validation cohorts, diagnostic labels were first mapped to the standardized Zhongshan ontology (detailed mapping rules for each external dataset are provided in [Tables S4–S11](#); notably, the Xiamen cohort shared an identical ontology with Zhongshan and required no mapping). Evaluation was strictly restricted to the intersection of label spaces shared

between the source domain and each respective target dataset (see [Table S12](#) for tasks performed on each dataset and [Tables S13-14](#) for detailed prevalence). Following this standardization process, the downstream diagnostic tasks were divided into the following three tiers: (1) ECG tasks: We benchmarked the model on 45 distinct diagnostic tasks organized into six major diagnostic categories: ischemia, arrhythmia, conduction abnormalities, structural changes, pacing status, and electrolytic disturbances. This evaluation spanned the internal Zhongshan dataset and seven external cohorts (Xiamen cohort and six public benchmarks), ensuring comprehensive validation across diverse populations and recording devices. (2) Echo tasks: To assess the model's capacity to infer cardiac structural and functional abnormalities solely from ECG waveforms, we defined 39 diagnostic tasks grouped into seven categories: normal function, myocardial, valvular, shunt, vessel, pericardial, and arrhythmia-related anomalies. Validation was conducted specifically on the Zhongshan, Shiyuan, and BIDMC cohorts, which provided the necessary paired echocardiographic ground truth. To ensure pathophysiological consistency between modalities, we restricted the analysis to ECG-echocardiography pairs recorded within a maximum interval of 30 days, with the time window typically ranging between 14 days and one month. (3) Clinical tasks: The identification of complex cardiomyopathies, specifically alcoholic cardiomyopathy (AC), noncompaction of ventricular myocardium (NCVM), Takotsubo cardiomyopathy (TC), arrhythmogenic right ventricular cardiomyopathy (ARVC), and cardiac amyloidosis (CA), was evaluated primarily on the Zhongshan cohort. In this setting, ground truth labels were rigorously established through multi-modal clinical synthesis (integrating imaging, laboratory biomarkers, and pathological biopsy) and verified by senior diagnosticians, ensuring a gold standard for these challenging conditions.

ECG-CLIP Framework and Training Objectives

We adopted the dual-objective framework established by the Merl architecture⁹, scaling it to a large-scale, longitudinal clinical cohort to derive robust and transferable representations directly from raw ECG signals and clinical narratives. The training process optimizes a composite loss function combining Cross-Modal Alignment (CMA, the left panel of [Figure 1B](#)) and Uni-Modal Alignment (UMA, the right

panel of Figure 1B).

Cross-Modal Alignment (CMA). This objective bridges the semantic gap between signal and text. Given a batch of N pairs consisting of an ECG signal (x) and its corresponding text report (t), distinct encoders, namely a ResNet-based ECG encoder (f_{ecg}) and a Med-CPT text encoder (f_{text}), project inputs into initial embeddings. These are subsequently mapped by non-linear projectors (g_{ecg} and g_{text}) into a shared latent space with normalized embeddings e_i and r_i . We employ a symmetric InfoNCE loss to maximize the cosine similarity between matched pairs while suppressing unmatched negatives:

$$L_{CMA} = -\frac{1}{2N} \sum_{i=1}^N \left(\log \frac{\exp(e_i \cdot r_i^T / \tau)}{\sum_{j=1}^N \exp(e_i \cdot r_j^T / \tau)} + \log \frac{\exp(r_i \cdot e_i^T / \tau)}{\sum_{j=1}^N \exp(r_i \cdot e_j^T / \tau)} \right) \quad (1)$$

where τ denotes the temperature hyper-parameter, set to 0.07. This mechanism ensures that the ECG representations capture the detailed diagnostic semantics inherent in the expert-written reports.

Uni-Modal Alignment (UMA). To further refine the robustness of the signal encoder without relying on input-level augmentations that risk distorting pathological waveforms, we employed the UMA strategy within the ECG domain. This approach utilizes dropout as a data augmentation technique. For a given ECG input, two independent dropout masks (with a ratio $p = 0.1$) are applied to the encoder's output, generating two distinct but semantically equivalent views (\hat{e}_i, \tilde{e}_i) of the same recording. A contrastive loss is then minimized to enforce consistency between these views:

$$L_{UMA} = -\frac{1}{2N} \sum_{i=1}^N \left(\log \frac{\exp(\hat{e}_i \cdot \tilde{e}_i^T / \tau)}{\sum_{j=1}^N \exp(\hat{e}_i \cdot \tilde{e}_j^T / \tau)} + \log \frac{\exp(\tilde{e}_i \cdot \hat{e}_i^T / \tau)}{\sum_{j=1}^N \exp(\tilde{e}_i \cdot \hat{e}_j^T / \tau)} \right) \quad (2)$$

By compelling the model to maintain representational invariance under latent perturbation, UMA encourages the learning of resilient, high-level features critical for downstream generalization.

Model Variants and Baseline Comparisons

While our framework builds upon the dual-objective architecture established by Merl, our study focuses on the effects of scaling this approach from the original approximately 0.8 million recordings (MIMIC-IV) to a large-scale cohort comprising approximately 2.8 million ECG-text pairs for pre-training.

Empirical results (Table S15) indicated that scaling to ResNet-50 yielded marginal gains (or even performance drop) compared to ResNet-34. Consequently, ECGCLIP-R34 was selected as our primary model. However, to ensure a fair comparison with the prior art, we explicitly included ECGCLIP-R18 and Random Init-R18 to benchmark directly against the Merl-R18 baseline, thereby isolating the benefits of data scaling from architectural depth.

Implementation Details

Signal Preprocessing. To ensure consistent model input across heterogeneous cohorts, we implemented a unified signal processing pipeline. The primary model input was defined as the eight independent ECG leads (I, II, V1–V6), as the remaining standard limb leads (III, aVR, aVL, aVF) are linear combinations of leads I and II. Raw signals underwent a third-order Butterworth band-pass filter (0.5–100 Hz) to mitigate baseline wander and high-frequency noise, followed by a 50 Hz notch filter (with harmonics) to eliminate power-line interference. Crucially, both filtering steps utilized a forward-backward approach to ensure zero-phase distortion, preserving the temporal fidelity of the original waveforms. Temporal and sampling consistency was strictly enforced. All recordings were resampled to a uniform 500 Hz. We standardized the input duration to a fixed 10-second window. For datasets with variable signal lengths, such as CPSC2018 (ranging from 6 to 60 seconds) or subsets of the Georgia cohort, recordings falling short of this duration were zero-padded, while those exceeding it were truncated to the initial 10-second segment. This resulted in a uniform (8, 5000) tensor for each preprocessed ECG signal. For baseline comparisons requiring standard 12-lead input (i.e., Merl), we reconstructed the four dependent leads (III, aVR, aVL, aVF) via Einthoven's law linear combinations when they were not explicitly available, ensuring fair comparison metrics.

Model Hyperparameters and Training Setup. The training pipeline comprised a pre-training phase followed by downstream fine-tuning. Pre-training was conducted on two NVIDIA A100-80GB GPUs using a mixed-precision framework for 20 epochs (pre-training hyperparameters are detailed in Table S16). We instantiated the ECG encoder with three variants of the ResNet-1D architecture (ResNet-18, -34, and -50). To maximize hardware utilization under memory constraints, per-GPU batch sizes were

set to 768, 700, and 512, respectively, resulting in effective global batch sizes of 1536, 1400, and 1024. To bridge the linguistic gap between the original Chinese diagnostic reports and the English-based Med-CPT text encoder, we utilized GPT-4o to translate all clinician-authored narratives into standardized English. The Med-CPT architecture was then initialized with these translated reports; to preserve learned semantic knowledge while adapting to ECG contexts, the initial nine layers were frozen, and only the final three layers were fine-tuned. Optimization was performed using AdamW with a Cosine Annealing with Warm Restarts scheduler. The restart period (T_0) was set to 5,000 iterations with a fixed cycle length ($T_{mult}=1$) and a minimum learning rate (η_{min}) of 1×10^{-8} . Checkpoints were selected based on the highest signal-report retrieval accuracy on the validation set. For downstream fine-tuning (see [Table S17](#) for detailed settings), the model underwent 4,000 iterations. In this phase, we employed a linear warmup strategy for the first 200 iterations (warmup to 1×10^{-6}), followed by standard cosine decay. Performance was monitored every 200 iterations, and the checkpoint achieving the highest macro-averaged Area Under the Precision-Recall Curve (PRAUC) on the validation set was selected for final evaluation on the independent test sets.

Interpretability and Feature Visualization

Saliency Mapping via Integrated Gradients. To identify specific ECG segments driving diagnostic predictions, we utilized the Integrated Gradients (IG) method, an axiomatic attribution technique that satisfies sensitivity and implementation invariance. We computed the path integral of the gradients for the predicted class score with respect to the input ECG signal, scaling from a zero-baseline reference to the actual input. To prioritize clinical relevance and visual clarity, we focused exclusively on positive attribution regions (features that positively contribute to the target diagnosis) by masking negative gradients ($IG < 0$). These rectified saliency maps were superimposed onto the original waveforms as heatmaps, highlighting the precise morphological patterns (e.g., ST-segment elevation, pathological Q waves) that the model identified as diagnostic evidence.

Manifold Visualization via t-SNE. To assess the global discriminative capability of the learned representations, we visualized the high-dimensional latent space using t-SNE. We extracted embeddings

from the pre-trained ECG encoder to evaluate the representation quality. To mitigate the visual ambiguity inherent in multi-label clinical cases, we restricted this analysis strictly to single-label samples (patients presenting with a solitary pathology), ensuring that cluster boundaries reflect distinct disease phenotypes rather than comorbidity overlap. We conducted a hierarchical visual analysis spanning from broad diagnostic categories, such as delineating electrolytic disturbances from pacing status, down to fine-grained intra-cluster distinctions, such as resolving JPB versus JEB. Stratified random sampling was employed to ensure balanced representation. For intra-cluster visualization, we applied a minimum count threshold ($N \geq 200$) to exclude sparse classes that might introduce noise, except rare disease categories, ensuring that low-prevalence cardiomyopathies remained visible. The t-SNE algorithm was configured with a perplexity of 30 and run for 1,000 iterations to guarantee stable convergence.

Evaluation Metrics and Statistical Analysis

Performance Metrics and Thresholding. Model performance was comprehensively evaluated using the Area Under the Precision-Recall Curve (PRAUC), Area Under the Receiver Operating Characteristic Curve (ROAUC), sensitivity, specificity, F1-score, and the Matthews Correlation Coefficient (MCC). PRAUC was reported as the primary endpoint given its robustness in evaluating multi-label classification tasks under significant class imbalance. For metrics requiring binary decisions (sensitivity, specificity, F1-score, MCC), operating thresholds were determined strictly on the internal validation set. Specifically, we identified the threshold that maximized the F1-score for each diagnostic task, and these fixed thresholds were subsequently applied uniformly to all external validation cohorts to ensure unbiased evaluation.

Statistical Significance and Uncertainty Quantification. To rigorously quantify statistical uncertainty, 95% Confidence Intervals (CI) were computed using nonparametric bootstrap resampling with 1,000 iterations. Random seeds were fixed to ensure reproducibility. Comparative analysis between models was conducted using a two-sided paired permutation test with 1,000 permutations. This non-parametric approach makes no assumptions about the underlying data distribution. In each permutation, the

predicted probabilities of the two comparison models were randomly swapped for each sample to generate a null distribution of performance differences. The P-value was calculated as the proportion of permuted differences exceeding the observed difference in absolute value. Statistical significance was defined as two-sided $P < 0.05$. All statistical analyses were performed using Python (v3.10) and scikit-learn (v1.5.0).

Code availability. The source code for model inference, pre-trained weights for ECGCLIP, and the scripts for reproducing the downstream evaluations are available on GitHub at <https://github.com/dt-yuhui/ECG-CLIP>.

Ethics and Data Security

The study was approved by the institutional review board (registration number B2024-497). All data collection and use adhered to institutional review board (IRB) protocols, with informed consent obtained for the use of patient data. Privacy-preserving techniques, including pseudonymization and secure data storage, were employed to protect patient confidentiality.

Results

Dataset composition and evaluation protocol

Our study utilized a multicenter dataset comprising 4,875,803 ECG studies from 2,094,169 patients (Figure 1A and Table S1-S2). The model architecture is shown in Figure 1B. To systematically quantify the clinical utility of ECGCLIP, we assessed its diagnostic performance across a comprehensive taxonomy comprising three distinct tiers: 45 standard ECG tasks, 39 ECHO tasks, and 5 comprehensive clinical tasks (Figures 1C–1D; Table S3).

Diagnostic performance across cardiac conditions diagnosed by ECG

Within the first tier of evaluation, 45 ECG tasks were clustered into six categories (ischemia, arrhythmia, conduction, structure, pacing, and electrolyte abnormality). While the primary Zhongshan cohort encompasses all 45 labels, evaluation on the seven external datasets was strictly conducted on their available subsets of these tasks (MIMIC-IV: 37; PTB-XL: 33; Georgia: 33; Chapman: 27; Xiamen: 26; UK Biobank: 24; CPSC: 8). Analysis of disease prevalence across cohorts revealed a distinct spectrum of disease frequency (Figure 2A, Tables S3, S13). Common arrhythmias, such as atrial fibrillation, atrial flutter, sinus bradycardia, and premature ventricular beats, predominate; conduction and structural heart diseases occur at intermediate rates; and pacing-related, electrolyte, and ischemic conditions are rare (<1%). This pronounced heterogeneity reflects both clinical diversity and class imbalance challenges in model training and evaluation.

We present a heatmap of PRAUC values (rows: labels, columns: cohorts) to visualize model robustness and generalizability, where uniformly high performance (indicated by deep red) denotes strong generalization, and color variation signals cohort-specific differences (Figure 2B; full results in Tables S18–S25). To quantify cross-population stability, we evaluated the average macro-performance of ECGCLIP against baseline models across all eight datasets (Figure 2C). Baseline models showed limited transferability, while ECGCLIP-R18 and ECGCLIP-R34 achieved much higher mean PRAUCs. Notably, ECGCLIP-R34 outperformed comparators across all validated cohorts except PTB-XL, demonstrating robust cross-dataset generalization. To further dissect the performance across these 45

tasks, we stratified the evaluation across their six predefined categories. By comparing the categorical average performance between the internal Zhongshan test set and a combined external test cohort, we demonstrated that ECGCLIP models yielded significant improvements across all categories except electrolyte abnormality (Figure 2D).

Having established the model's broad generalizability, we next sought to quantify the specific diagnostic gains attributable to our pre-training strategy. We highlighted the top 22 ECG tasks where ECGCLIP-R34 demonstrated the most substantial PRAUC improvements relative to the Random Init-R18 baseline on Zhongshan test set (Figure 3). ECGCLIP-based models achieved consistently higher PRAUC values than both random initialization (Random Init-R18) and the Merl-R18 baseline. Among all architectures, ECGCLIP-R34 delivered optimal performance and exhibited clear improvements compared to the two control methods. For common and clinically significant arrhythmias, including AF (Atrial Fibrillation), VPB (Ventricular Premature Beats), AFL (Atrial Flutter), and SVT (Supraventricular Tachycardia), ECGCLIP-R34 maintained stable and superior PRAUC values relative to Random Init-R18 and Merl-R18, confirming its reliable performance in standard clinical ECG interpretation tasks. In AF, ECGCLIP-R34 achieved a PRAUC of 0.900 (95% CI: 0.892–0.907), outperforming Random Init-R18 (0.830) and Merl-R18 (0.839). In VPB, ECGCLIP-R34 reached a PRAUC of 0.913 (95% CI: 0.907–0.919), exceeding Random Init-R18 (0.768) and Merl-R18 (0.795). In AFL, ECGCLIP-R34 delivered a PRAUC of 0.821 (95% CI: 0.806–0.836), maintaining a lead over Random Init-R18 (0.739) and Merl-R18 (0.767). In SVT, ECGCLIP-R34 achieved a PRAUC of 0.757 (95% CI: 0.682–0.822), outperforming Random Init-R18 (0.631) and Merl-R18 (0.651). In particular, for the high-stakes task of STEMI, the performance gain of ECGCLIP-R34 was pronounced to achieve a PRAUC of 0.383 (95% CI: 0.332–0.436), representing a near-doubling of performance compared to both Random Init-R18 (0.199, 95% CI: 0.164–0.240) and Merl-R18 (0.206, 95% CI: 0.169–0.248).

Notably, the most substantial relative improvements in PRAUC were observed in low-prevalence or diagnostically challenging arrhythmias and conduction disorders, such as AAI, 2nd/3rd degree AVB, ATach, and VEB. In AAI, ECGCLIP-R34 achieved a PRAUC of 0.705 (95% CI: 0.639–0.764), representing a substantial increase compared with Random Init-R18 (0.190) and Merl-R18 (0.360). In

ATach, ECGCLIP-R34 delivered a PRAUC of 0.222 (95% CI: 0.194–0.252), while Random Init-R18 (0.093) and Merl-R18 (0.108) showed limited discriminative ability. In VEB, ECGCLIP-R34 achieved a PRAUC of 0.201 (95% CI: 0.138–0.266), outperforming Random Init-R18 (0.100) and Merl-R18 (0.092). For the remaining conditions, including APB, DDD, VAT, JEB, JTach, VVI, RVH, LAFB, OMI, Type II second degree AVB, QT Prolong, and IVB, ECGCLIP-R34 also achieved consistently higher PRAUC values than both Random Init-R18 and Merl-R18, further validating the generalized advantage of the ECGCLIP pre-training framework.

Diagnostic performance across cardiac conditions diagnosed by ECHO

We validated ECGCLIP on three cohorts: Zhongshan (39 labels), Shiyuan (36 labels), and BIDMC (17 labels) for ECHO-related tasks. Task labels were grouped into seven categories (normal, myocardial, valvular, shunt, vessel, pericardial, and arrhythmia abnormality) with varying prevalence (Figure 4A; Tables S3, S14). A comprehensive heatmap of PRAUC values demonstrated ECGCLIP’s enhanced overall diagnostic efficacy (Figure 4B; Tables S18, S26–S27). To quantify cross-cohort stability, we evaluated macro-performance across the three datasets (Figure 4C). A consistent performance hierarchy emerged (ECGCLIP-R34 > ECGCLIP-R18 > Merl-R18 > Random Init-R18), validating ECGCLIP’s robust generalization (Figure 4C). By further dissecting performance across the seven categories in the Zhongshan test set and combined external test set (Figure 4D), a similar trend was observed that deeper ECGCLIP networks consistently surpassed baselines across all categories.

Similar to the ECG tasks, we next sought to quantify the specific diagnostic gains attributable to our pre-training strategy for ECHO tasks within the Zhongshan test set (Figure 5). Across most diagnostic tasks, models initialized with the ECGCLIP framework consistently outperformed those with random initialization (Random Init-R18) and the Merl-R18. Notably, the most substantial relative improvements in PRAUC were observed in low-prevalence diseases that are not reliably diagnosed by conventional ECG, particularly EA, constrictive pericarditis (CP), and Dex. In EA, ECGCLIP-R34 achieved a PRAUC of 0.253 (95% CI: 0.089–0.433), representing a substantial relative increase compared with Random Init-R18 (0.003) and Merl-R18 (0.004), both of which performed near chance

level. In CP, ECGCLIP-R34 yielded a PRAUC of 0.175 (95% CI: 0.096–0.283), substantially outperforming Random Init-R18 (0.043) and Merl-R18 (0.037). In Dex, ECGCLIP-R34 reached a PRAUC of 0.121 (95% CI: 0.004–0.312), while Random Init-R18 (0.000) and Merl-R18 (0.004) showed negligible discriminative ability. For the remaining conditions, including oHCM, CHD, VA, VSD, BAE, RMWSSFLV, AD, and AVC, ECGCLIP-R34 also achieved consistently higher PRAUC values than both Random Init-R18 and Merl-R18, further validating the generalized advantage of the ECGCLIP pre-training framework.

Training efficiency and interpretability

Evaluating the average performance across all eight datasets (Figure 6A), ECGCLIP models consistently outperformed alternative approaches across all data regimes (1%, 10%, and 100% training data fractions). Notably, ECGCLIP-R34 achieved the highest overall PRAUC. This superiority held true across individual cohorts with the sole exception of the PTB-XL dataset at higher data fractions (10% and 100%) (Figure S1; Table S28). Notably, ECGCLIP-R34 fine-tuned on just 10% of the data matched or exceeded the performance of baseline models trained on the entire 100% dataset, showing strong sample efficiency and transfer learning. These performance gaps persisted as data increased to 100%. Conversely, the randomly initialized baseline (Random Init-R18) exhibited the lowest performance using 1% training data, exposing limitations of standard supervised training from scratch in data-scarce scenarios (Figure S1).

To understand the underlying representational quality driving this data efficiency, we visualized the learned embeddings using t-SNE. ECGCLIP-R34 yielded more distinct and compact t-SNE clusters than Merl. This structural delineation was evident both at the macro-level, where the six overarching clinical categories were clearly delineated (Figure 6B), and at the micro-level, successfully separating specific phenotypes within complex domains such as arrhythmias (Figure 6C) and conduction abnormalities (Figure 6D). In contrast, the Merl representations displayed heavy inter-class overlap. These results confirm that ECGCLIP pre-training learns clinically meaningful, well-separated representations (Figures S2–S3).

Having established the global validity of the latent space via t-SNE, we next evaluated spatial transparency at the waveform level using Integrated Gradients (IG). These saliency heatmaps successfully localized the physiological basis of model predictions. For STEMI, the model focuses on diagnostic ST–T segments and elevation patterns (Figure 7A). For APC, it attends to premature P-waves and associated QRS-T complexes (Figure 7B). For atrial fibrillation, it captures global rhythm disturbances consistent with clinical criteria (Figure 7C). For AFL, it identifies characteristic sawtooth flutter waves (Figure 7D). For LBBB and RBBB, attention aligns with their defining QRS morphologies (Figure 7E–F). For DDD pacing, it focuses on pacing spikes and abnormal QRS complexes (Figure 7G). For pre-excitation, it highlights delta waves in the early QRS complex (Figure 7H). Interpretability results across cardiac conditions are summarized in Figures S4–S16.

Potential of recognizing ECG characteristics on rare heart disease

The t-SNE visualization shows ECGCLIP-R34 forms more distinct, compact clusters for rare diseases than the Merl model, with clearer class boundaries and stronger discriminative features, especially for cardiac amyloidosis (CA) (Figure 8A). For CA detection, PRAUC values were: Random Init-R18 0.047, Merl-R18 0.049, ECGCLIP-R18 0.165, ECGCLIP-R34 0.201 (Figure 8B). ECGCLIP models yielded meaningful improvements (Figure 8C), with ECGCLIP-R18 achieving approximately 5-fold higher PRAUC over baselines, which performed near chance. Beyond discriminative accuracy, attention heatmaps confirmed the model focuses on clinically meaningful ECG regions (Figure 8C). Saliency maps consistently highlight P-wave and T-wave segments, which are established electrocardiographic markers of diastolic function. This physiological alignment is particularly noteworthy given that diastolic dysfunction is a hallmark clinical manifestation of CA^{10,11}.

Ablation analysis of the multimodal pre-training framework

To systematically isolate the architectural and data-driven components contributing to ECGCLIP’s diagnostic performance, we conducted a comprehensive ablation analysis evaluating both pre-training scale and model capacity across the three diagnostic tiers (Table S15).

First, to quantify the impact of our large-scale, high-fidelity supervision, we compared ECGCLIP-R18 against the baseline Merl-R18. Both models share an identical ResNet-18 contrastive learning backbone. However, Merl was pre-trained on MIMIC-IV-ECG (around 800k pairs), whereas ECGCLIP-R18 leveraged our extensive cohort of over 2.8 million expert-annotated pairs. Across all evaluations, ECGCLIP-R18 consistently outperformed Merl-R18. For standard ECG interpretation and ECHO diagnosis, ECGCLIP-R18 demonstrated robust baseline improvements (average PRAUC 0.388 vs. 0.370, and 0.283 vs. 0.273, respectively). Crucially, the advantage of this high-fidelity supervision was most evident in the most challenging rare disease diagnosis. While absolute performance metrics for these ultra-rare conditions remain inherently modest, ECGCLIP-R18 achieved a PRAUC of 0.035, more than tripling the performance of Merl-R18 (0.011), which performed near random chance.

Second, we investigated the effect of network depth on representational power by scaling the ECGCLIP backbone across ResNet-18, ResNet-34, and ResNet-50 architectures. Increasing the model capacity from ResNet-18 to ResNet-34 yielded systematic PRAUC improvements across all diagnostic tiers. The average PRAUC increased from 0.388 to 0.409 in ECG interpretation, from 0.283 to 0.292 in ECHO diagnosis, and reached a peak of 0.045 in rare disease diagnosis. However, further scaling the architecture to ResNet-50 did not yield additional benefits. In fact, performance degraded across multiple tasks, most notably in rare disease diagnosis where the PRAUC dropped to 0.028.

Discussion

To our knowledge, ECGCLIP represents a substantial advance over prior AI-ECG systems, driven primarily by leveraging large-scale data, clinical diversity, and high-fidelity multimodal alignment. Trained on 2,781,289 expert-authored ECG–text pairs from 1,298,359 patients, it enables robust, general-purpose cardiac assessment from a single standard ECG, a capability that smaller or task-specific models typically struggle to achieve. We validate its demographic robustness and broad generalizability across nine external cohorts spanning distinct geographies and care contexts. Notably, ECGCLIP excels not only on standard diagnostic tasks but also on challenging clinical targets, such as screening for occult structural phenotypes (e.g., HFrEF), underscoring its potential to streamline workflows and augment clinical decision-making. These analyses demonstrate that ECGCLIP learns clinically meaningful global representations, and its local feature attribution is physiologically coherent with established electrocardiographic criteria, mitigating the “black box” concern.

Contextualizing Within Existing Research

The evolution of AI-enabled electrocardiography began with task-specific models optimized for single clinical targets. Early algorithmic frameworks, such as ECG-ResNet⁵, achieved leading performance for specific rhythm classifications but lacked comprehensive clinical deployment. Subsequently, systems like Cardiologs¹⁴ provided high-level clinical evidence for single-disease alerts, demonstrating through randomized controlled trials (RCTs) that AI-detected atrial fibrillation translates into actionable clinical interventions (sensitivity 97–98%, specificity 98–99%). Similarly, groundbreaking single-task models established the potential of AI-ECG² to detect asymptomatic left ventricular dysfunction with an AUROC of 0.932. Extending this utility to structural valvular diseases, Cohen-Shelly et al.¹⁶ developed a convolutional neural network for AI-enabled ECG screening of aortic stenosis (AS), achieving an AUC of 0.85 in a large cohort. This work demonstrated that AI-ECG could effectively identify structural valvular heart disease, serving as a powerful community screening tool and highlighting the clinical utility of single-task models for actionable cardiovascular conditions. More recently, Liang et al.¹⁵ extended this focus to regurgitant valvular heart diseases to predict future mitral and tricuspid

regurgitation with robust prognostic performance (C-indices of 0.774 and 0.793, respectively). While these pioneering models proved that discriminative accuracy could drive meaningful clinical action, they inherently function as isolated, narrow classifiers, lacking the flexibility to serve as panoramic screening tools for the full spectrum of cardiovascular disease.

To overcome these single-disease silos, recent landmark efforts have established the feasibility of foundation models in cardiovascular medicine. Zhou et al.¹² first demonstrated the viability of zero-shot ECG diagnosis using natural language supervision for 18 common cardiac conditions. Subsequently, Tian et al.⁷ proposed KED which advanced this paradigm by incorporating medical knowledge enhancement, achieving robust zero-shot generalization across diverse populations for standard rhythm and morphological abnormalities. Concurrently, DeepECG⁸ leveraged automated free-text label extraction on a massive scale, demonstrating high-throughput scalability across 77 conventional ECG tasks. These foundational models validated the power of large-scale pre-training; however, they primarily focused on common electrophysiological conditions and relied heavily on automatically generated or ICD-derived labels.

Building upon this landscape, the superiority of ECGCLIP over prior comparative models can be conceptualized across three key dimensions: supervision quality, diagnostic breadth, and clinical depth. First, regarding supervision quality, unlike prior systems dependent on automated label extraction, ECGCLIP was trained exclusively on approximately 3 million clinician-authored interpretations from a specialized ECG laboratory. This semantic richness captures nuanced language that encodes contextual reasoning (e.g., “new ST concordance in LBBB”) and diagnostic uncertainty (e.g., “possible anterior ischemia”). Second, regarding diagnostic breadth, ECGCLIP breaks single-disease silos to provide a unified assessment across 89 distinct tasks. For instance, while a prior single-task model achieved an AUROC of 0.932 for detecting low ejection fraction², our unified framework yields comparable discriminative precision (AUROC 0.96 for HFrEF [EF <40%]) alongside 88 other diagnoses. Third, the model successfully extends the predictive boundaries of the standard ECG to include complex echocardiography-derived structural phenotypes and elusive rare cardiomyopathies. In the detection of cardiac amyloidosis, for example, model performance improved significantly from near-random chance

(PRAUC around 0.05) to 0.201. By addressing these clinically challenging targets through expert-level multimodal alignment, ECGCLIP establishes a new paradigm for panoramic cardiovascular screening.

Opportunistic Screening for Structural and Rare Phenotypes

The standard ECG is inexpensive, ubiquitous, and routinely acquired, yet it is traditionally considered “blind” to structural heart disease. Conversely, echocardiography is information-rich but costly, scarce, and often delayed. By enabling the robust prediction of echo-level phenotypes, including HFrEF, valvular lesions (AS, MR), and even rare conditions like cardiac amyloidosis directly from ECG waveforms, ECGCLIP transforms routine ECG into an opportunistic structural heart screen^{2,15,16}. In this context, the model acts as a highly scalable gatekeeper for primary care. Rather than functioning as a standalone diagnostic tool, it effectively enriches the pretest probability of high-risk patients, thereby optimizing the allocation of scarce and expensive confirmatory resources such as echocardiography. This capability has profound implications for early detection, especially in resource-limited settings where echocardiography access is limited but ECGs are abundant. Notably, ECGCLIP-R34 demonstrated robust performance in both common and challenging diagnostic scenarios. For routine cardiovascular conditions such as ASD, HFrEF, and TR, all models achieved moderate-to-high PRAUC values, with ECGCLIP-R34 consistently delivering the best results. However, the most notable relative advantages of ECGCLIP-R34 were observed in diseases with low prevalence and conditions that are not routinely diagnosed via standard ECG, such as EA, CP, and Dex. While baselines performed near chance levels (PRAUC < 0.05), ECGCLIP-R34 achieved multi-fold PRAUC improvements. This suggests that the rich, task-agnostic representations learned through multimodal pre-training enable the model to capture subtle, non-obvious ECG patterns that are critical for identifying these rare and diagnostically elusive conditions. The clinical value is particularly evident in challenging, high-stakes diagnoses like cardiac amyloidosis, a condition often missed until advanced stages. Here, the proposed framework achieves a three-fold improvement in PRAUC over baseline models. Furthermore, saliency maps highlight diastolic-relevant P-wave and T-wave regions, consistent with the known pathophysiology of the disease. Similarly, constrictive pericarditis and Ebstein anomaly, rare and

electrocardiographically subtle, show substantial relative gains, underscoring ECGCLIP's ability to extract latent signals invisible to both clinicians and prior AI systems¹⁷.

Drivers of Data Efficiency, Global Generalizability, and Health Equity

The robust diagnostic capabilities and high data efficiency of ECGCLIP are fundamentally driven by our multimodal pre-training strategy. As demonstrated in the ablation analysis, utilizing large-scale, clinically rich text supervision is the primary catalyst for capturing elusive cardiovascular conditions, significantly outperforming the Merl-R18 baseline trained on smaller corpora. Additionally, optimizing the network capacity to a ResNet-34 architecture provides an ideal balance between representational depth and generalization, effectively preventing the overfitting and representational saturation observed in deeper architectures like ResNet-50. Because this optimized foundation extracts highly stable and transferable cardiac features, it substantially reduces the demand for extensive downstream fine-tuning. Consequently, ECGCLIP matches or exceeds the performance of baselines trained on the entire dataset using only a 10% fraction of the downstream training data, and it remains effective even in the 1% data regime (Figure 6A, Figure S1, Table S28).

Driven by these universally learned representations, ECGCLIP demonstrates robust cross-population generalizability. Across nine independent external cohorts spanning North America, Europe, and Asia (Table S1), the model consistently outperformed baselines. Notably, it maintained high diagnostic accuracy across 45 ECG tasks and 39 echocardiographic phenotypes in non-Asian demographics (Table S15). Unlike traditional algorithms that frequently overfit to the demographic distribution of their primary training data, ECGCLIP captures physiologically meaningful and population-agnostic features of cardiac dysfunction. Ultimately, this dual advantage of high data efficiency and global generalizability actively promotes health equity. Because the model requires minimal local data to achieve high diagnostic accuracy and performs equitably across diverse international cohorts, small clinics or low-resource regions can easily adapt the pre-trained framework to their specific populations. This helps mitigate potential hidden demographic biases and democratizes access to advanced cardiovascular diagnostics for historically marginalized groups.

Limitations

Despite the robust diagnostic performance of ECGCLIP, several limitations warrant consideration. First, the large pre-training dataset is predominantly sourced from a single institution comprising a Chinese population. Although we successfully validated the model across diverse international cohorts spanning Asia, North America, and Europe, the geographical concentration of the training data necessitates cautious extrapolation. The model's generalizability to entirely unrepresented demographics, such as African and Latin American populations who may exhibit distinct baseline electrocardiographic phenotypes, requires further dedicated evaluation. Furthermore, extreme age groups including pediatric and geriatric patients remain underrepresented. Second, the evaluation of echocardiographic and rare disease tasks was strictly confined to cohorts possessing paired clinical gold standards like echocardiograms or pathological biopsies. Because such comprehensive ground-truth data are scarce in routine primary care, the model's generalizability to less extensively phenotyped community populations remains unverified. Third, although ECGCLIP achieved substantial improvements in detecting rare diseases compared to baselines, absolute PRAUC values remained inherently modest. In a real-world primary care setting, this could translate to a high rate of false positives, potentially leading to alert fatigue and unwarranted downstream testing. Future iterations must optimize the precision-recall trade-off to ensure cost-effectiveness. Fourth, ECGCLIP performed on par with a fully supervised baseline on the PTB-XL dataset, exposing a limitation of our multimodal pre-training paradigm. By biasing representations toward complex morphological pathologies (Table S23, such as 3° AVB and RAH), the framework may compromise sensitivity to simple deterministic rules and hardware artifacts. Consequently, ECGCLIP underperformed the baseline in detecting artificial pacing spikes (e.g., VP) and basic rate classifications (e.g., SBrad). Future iterations should aim to address this bias to ensure uniform sensitivity across both high-level clinical phenotypes and low-level electrophysiological signals. Finally, the retrospective design of our study mandates careful examination regarding clinical workflow integration. Prospective randomized controlled trials are imperative to determine whether the improved discriminative accuracy of ECGCLIP translates into meaningful clinical outcomes, such as

reduced time-to-diagnosis or lower cardiovascular mortality.

Conclusion

ECGCLIP is a meaningful step forward in cardiac AI, trained on the largest expert-curated multimodal ECG dataset to date, achieving robust generalizability across 89 tasks. Its key strength is excelling at routine diagnostics while improving performance on low-prevalence, challenging conditions, transitioning from single-disease classifiers to a unified panoramic screening system for comprehensive cardiovascular assessment. With high data efficiency, interpretability and cross-cohort generalization, it redefines the ECG as a gateway to holistic cardiac health, with limitations in population diversity and rare disease samples that enable iterative improvement for clinical integration and global diagnostic democratization.

Acknowledgments

The authors would like to acknowledge the expert technical assistance of the data center staff at Zhongshan Hospital, Fudan University.

Sources of Funding

This work was supported by the National Natural Science Foundation of China (Grant G2021-046), Shanghai Pujiang Talent Program (24PJD014), and Fudan AI4S program (FudanX24A1056).

Disclosures

Dr. Hongyang Lu is an employee of Medtronic; however, Medtronic did not sponsor this study. All other authors declared no conflicts of interest.

Reference

1. Kligfield P, Gettes LS, Bailey JJ, et al. Recommendations for the standardization and interpretation of the electrocardiogram: part I: The electrocardiogram and its technology: a scientific statement from the American Heart Association Electrocardiography and Arrhythmias Committee, Council on Clinical Cardiology; the American College of Cardiology Foundation; and the Heart Rhythm Society: endorsed by the International Society for Computerized Electrocardiology. *Circulation*. 2007;115(10):1306-1324. doi:10.1161/CIRCULATIONAHA.106.180200
2. Attia ZI, Kapa S, Lopez-Jimenez F, et al. Screening for cardiac contractile dysfunction using an artificial intelligence-enabled electrocardiogram. *Nat Med*. 2019;25(1):70-74. doi:10.1038/s41591-018-0240-2
3. Siontis KC, Noseworthy PA, Attia ZI, Friedman PA. Artificial intelligence-enhanced electrocardiography in cardiovascular disease management. *Nat Rev Cardiol*. 2021;18(7):465-478. doi:10.1038/s41569-020-00503-2
4. Attia ZI, Harmon DM, Behr ER, Friedman PA. Application of artificial intelligence to the electrocardiogram. *Eur Heart J*. 2021;42(46):4717-4730. doi:10.1093/eurheartj/ehab649
5. Pham H, Egorov K, Kazakov A, Budenny S. Machine learning-based detection of cardiovascular disease using ECG signals: performance vs. complexity. *Front Cardiovasc Med*. 2023;10:1229743. doi:10.3389/fcvm.2023.1229743
6. Christopoulos G, Graff-Radford J, Lopez CL, et al. Artificial Intelligence-Electrocardiography to Predict Incident Atrial Fibrillation: A Population-Based Study. *Circ Arrhythm Electrophysiol*. 2020;13(12):e009355. doi:10.1161/CIRCEP.120.009355
7. Tian Y, Li Z, Jin Y, et al. Foundation model of ECG diagnosis: Diagnostics and explanations of any form and rhythm on ECG. *Cell Rep Med*. 2024;5(12):101875. doi:10.1016/j.xcrm.2024.101875
8. Nolin-Lapalme A, Sowa A, Delfrate J, et al. Foundation models for electrocardiogram interpretation: clinical implications. *Eur Heart J*. Published online January 22, 2026:ehaf1119. doi:10.1093/eurheartj/ehaf1119
9. Liu C, Wan Z, Ouyang C, Shah A, Bai W, Arcucci R. Zero-Shot ECG Classification with Multimodal Learning and Test-time Clinical Knowledge Enhancement. Published online 2024.
10. Liu D, Hu K, Niemann M, et al. Effect of combined systolic and diastolic functional parameter assessment for differentiation of cardiac amyloidosis from other causes of concentric left ventricular hypertrophy. *Circ Cardiovasc Imaging*. 2013;6(6):1066-1072. doi:10.1161/CIRCIMAGING.113.000683
11. Rapezzi C, Aimo A, Barison A, et al. Restrictive cardiomyopathy: definition and diagnosis. *Eur Heart J*. 2022;43(45):4679-4693. doi:10.1093/eurheartj/ehac543
12. Zhou X, Li T, Hayama H, et al. Diagnosis of cardiac conditions from 12-lead electrocardiogram through natural language supervision. *NPJ Digit Med*. 2025;8(1):697. doi:10.1038/s41746-025-02074-3
13. Yao X, Rushlow DR, Inselman JW, et al. Artificial intelligence-enabled electrocardiograms for identification of patients with low ejection fraction: a pragmatic, randomized clinical trial. *Nat Med*. 2021;27(5):815-819. doi:10.1038/s41591-021-01335-4

14. Liu WT, Lin C, Lee CC, et al. Artificial Intelligence-Enabled ECGs for Atrial Fibrillation Identification and Enhanced Oral Anticoagulant Adoption: A Pragmatic Randomized Clinical Trial. *J Am Heart Assoc.* 2025;14(14):e042106. doi:10.1161/JAHA.125.042106
15. Liang Y, Sau A, Zeidaabadi B, et al. Artificial intelligence-enhanced electrocardiography to predict regurgitant valvular heart diseases: an international study. *Eur Heart J.* 2025;46(44):4823-4837. doi:10.1093/eurheartj/ehaf448
16. Cohen-Shelly M, Attia ZI, Friedman PA, et al. Electrocardiogram screening for aortic valve stenosis using artificial intelligence. *Eur Heart J.* 2021;42(30):2885-2896. doi:10.1093/eurheartj/ehab153
17. Schrutka L, Anner P, Agibetov A, et al. Machine learning-derived electrocardiographic algorithm for the detection of cardiac amyloidosis. *Heart.* 2022;108(14):1137-1147. doi:10.1136/heartjnl-2021-319846

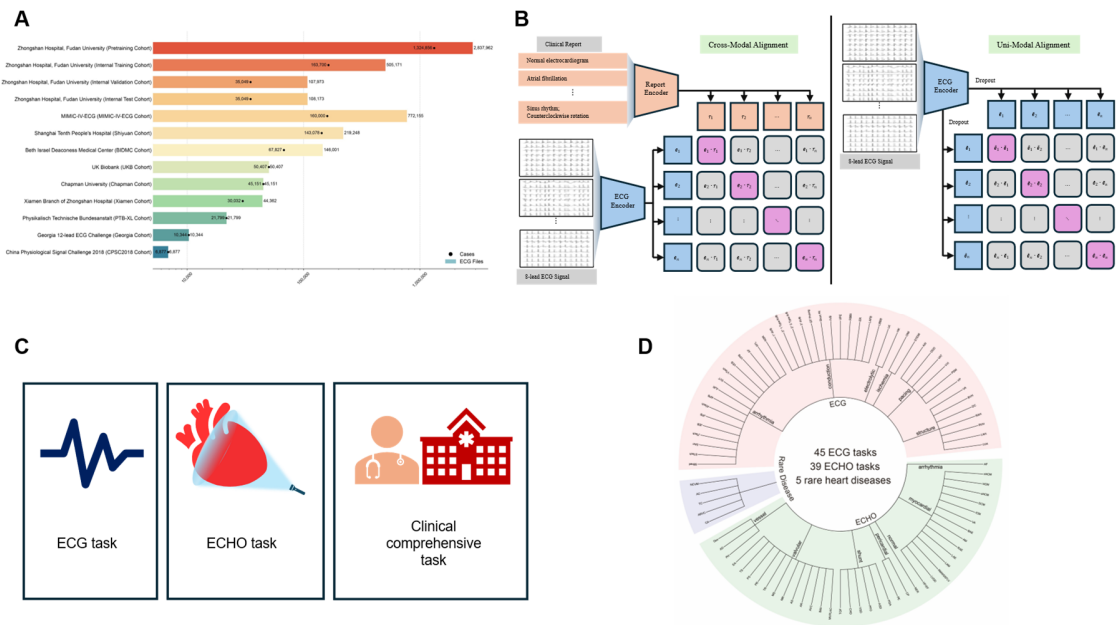


Figure 1. Data overview, model framework and downstream task taxonomy

A: Composition and scale of electrocardiogram (ECG) datasets across cohorts, summarizing the number of ECG files (bars) and unique clinical cases (points) for each cohort used in the study. The Zhongshan Hospital pretraining cohort provided the largest source of data, from which internal training, validation, and test splits were derived. Nine independent external cohorts were utilized to comprehensively evaluate model generalizability. The diverse, multi-source framework ensures robust development and validation across different patient populations and clinical settings.

B: ECGCLIP pretraining framework with dual alignment objectives. Cross-modal alignment (CMA, left) employs a report encoder to map clinical text (e.g., "Atrial fibrillation") into an embedding space where dot products between report representations ($R_1 \dots R_n$) and corresponding ECG representations ($E_1 \dots E_n$) from the ECG encoder are maximized. Uni-modal alignment (UMA, right) further enhances the robustness of ECG representations by utilizing dropout to generate two distinct but semantically equivalent views of the same recording, maximizing their similarity within each training batch.

C: Downstream tasks are stratified into three tiers based on diagnostic complexity. These include standard ECG tasks, echocardiographic tasks, and comprehensive clinical tasks targeting rare heart diseases.

D: Cardiovascular task taxonomy. The framework encompasses 45 ECG diagnostic tasks (categorized into ischemia, arrhythmia, conduction, structure, pacing, and electrolyte abnormality) and 39 echocardiogram (ECHO) tasks (covering normal function, myocardial, valvular, shunt, vessel, pericardial, and arrhythmia-related anomalies), alongside dedicated coverage of 5 rare cardiac diseases. This comprehensive taxonomy enables a panoramic evaluation across electrical and structural heart pathologies.

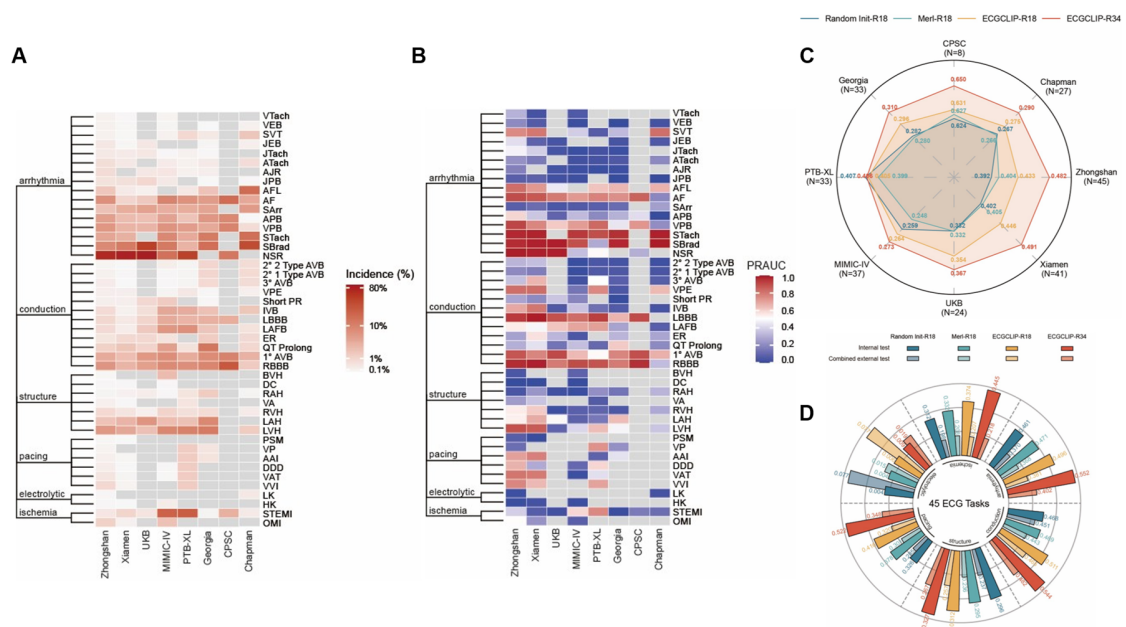


Figure 2. Cross-cohort generalizability and diagnostic performance of ECGCLIP across 45 ECG tasks.

A, Incidence heatmap. The heatmap displays the prevalence (percentages) of 45 distinct ECG diagnostic tasks, grouped into six categories (arrhythmia, conduction, structure, pacing, electrolytic, and ischemia), across eight independent cohorts (Zhongshan, Xiamen, UKB, MIMIC-IV, PTB-XL, Georgia, CPSC, and Chapman). Color intensity corresponds to incidence, with deep red indicating high prevalence and white indicating near-zero prevalence. Tasks unavailable in specific external cohorts are denoted by gray cells (number of evaluated tasks: Zhongshan, 45; MIMIC-IV, 37; PTB-XL, 33; Georgia, 33; Chapman, 27; Xiamen, 26; UKB, 24; CPSC, 8), highlighting the significant class imbalance and label heterogeneity across datasets.

B, PRAUC performance heatmap. This heatmap presents the PRAUC values for the ECGCLIP-R34 model across the same 45 tasks and eight cohorts. Color intensity represents PRAUC, with deep red indicating higher performance and blue indicating lower precision, and gray indicating unsupported tasks. While ECGCLIP-R34 achieves robust, high PRAUC values (warm colors) for common cardiovascular conditions, the absolute PRAUC values for low-prevalence tasks naturally remain more modest (cooler colors). This visual distribution accurately reflects the extreme class imbalance and the inherent clinical difficulty of screening for these rare targets.

C, Radar plot of cross-cohort performance. The radar plot compares the mean PRAUC of four models (Random Init-R18, Merl-R18, ECGCLIP-R18, ECGCLIP-R34) across eight external validation cohorts. The radius of each polygon represents the mean PRAUC for that model. ECGCLIP models, particularly ECGCLIP-R34, achieve larger polygon areas, demonstrating robust cross-cohort generalization.

D, Circular bar plot of internal versus external performance. This circular bar chart compares the PRAUC values of all four models on the internal test set (darker bars) and the combined external test set (lighter bars, generated by merging all seven external ECG cohorts) across 45 ECG tasks. Tasks are grouped by diagnostic category. For most tasks, ECGCLIP models maintain high and consistent PRAUC

values between internal and external testing, confirming robust generalization to unseen patient populations.

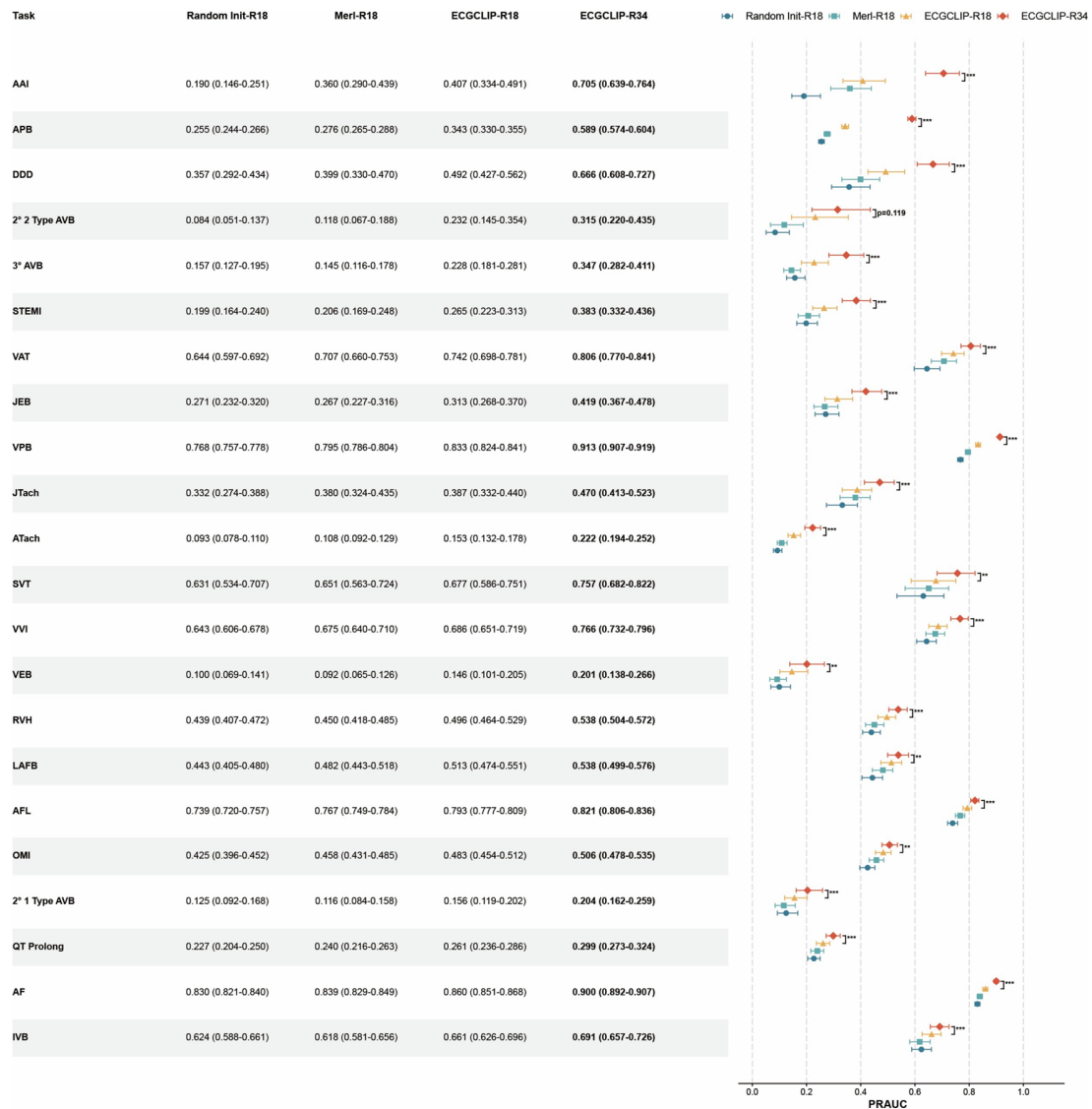


Figure 3. Performance comparison on the top 22 ECG diagnostic tasks showing the greatest improvement.

This figure presents a combined data table (left) and forest plot (right) illustrating model performance across 22 specific electrocardiographic tasks. These tasks were selected based on the greatest absolute PRAUC improvements achieved by ECGCLIP-R34 compared to the Random Init-R18 baseline. The left panel details the exact mean PRAUC alongside 95% confidence intervals (inside parentheses) for four model variants: Random Init-R18, Merl-R18, ECGCLIP-R18, and ECGCLIP-R34. The corresponding right panel visualizes these metrics, where points represent the mean PRAUC and error bars indicate the 95% confidence intervals derived from 1,000 bootstrap resamples. ECGCLIP models demonstrate substantial diagnostic improvements, particularly for low-prevalence conditions (for example, improving DDD from 0.357 to 0.666, and AAI from 0.190 to 0.705). Statistical significance between ECGCLIP-R34 and the second-best performing model for each task is indicated adjacent to the error bars by exact P values or asterisks (* $P < 0.05$, ** $P < 0.01$, *** $P < 0.001$).

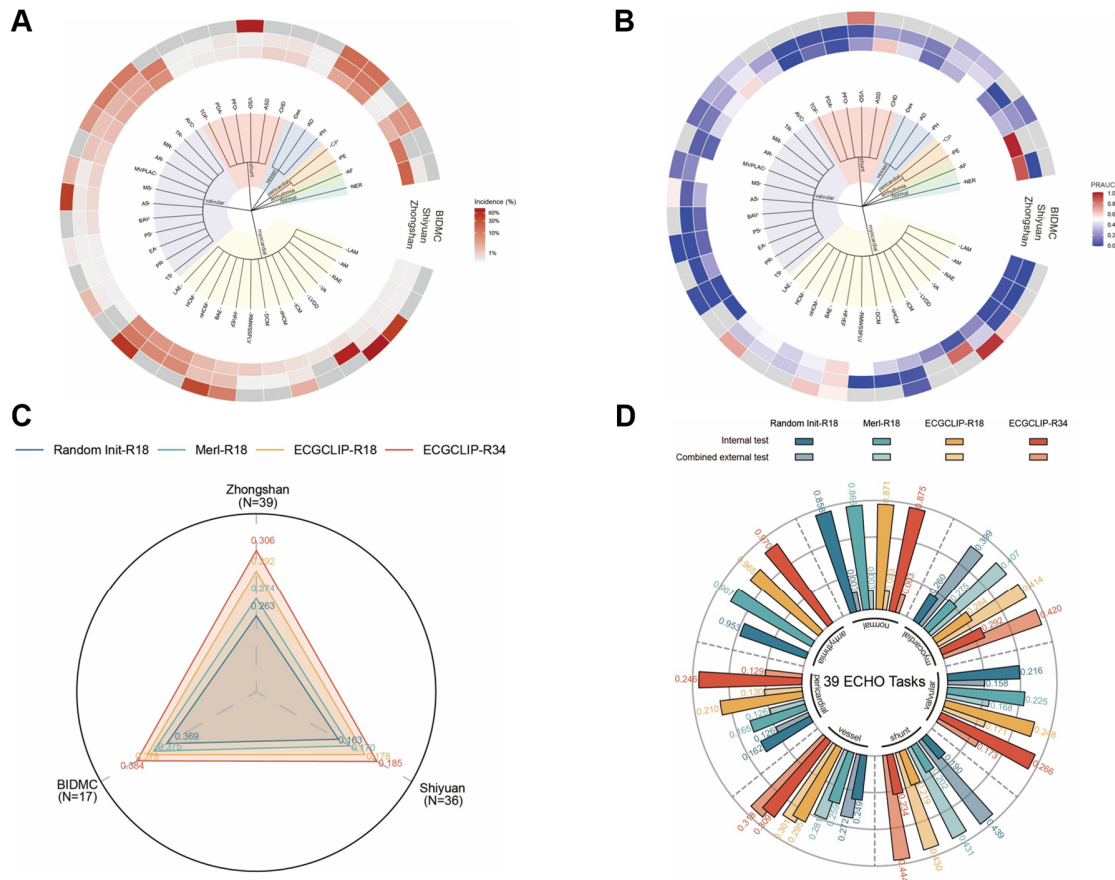


Figure 4. Diagnostic performance and cross-cohort generalizability of ECGCLIP across 39 ECHO tasks.

A: Incidence heatmap. The circular heatmap displays the prevalence (percentages) of 39 distinct echocardiographic diagnostic tasks, grouped into seven categories (myocardial, valvular, shunt, vessel, pericardial, normal function, and arrhythmia related anomalies), across three independent cohorts (Zhongshan, Shiyuan, and BIDMC). Color intensity corresponds to incidence, with deep red indicating high prevalence and white indicating near-zero prevalence. Tasks unavailable in specific external cohorts are denoted by gray cells (number of evaluated tasks: Zhongshan, 39; Shiyuan, 36; BIDMC, 17), highlighting the significant class imbalance across tasks and cohorts.

B: PRAUC performance heatmap. This circular heatmap presents the PRAUC values for the ECGCLIP-R34 model across the 39 tasks and three cohorts. Color intensity represents PRAUC, with deep red indicating higher performance, blue indicating lower precision, and gray indicating unsupported tasks. While ECGCLIP-R34 achieves robust, high PRAUC values (warm colors) for common structural phenotypes, the absolute PRAUC values for low-prevalence tasks naturally remain more modest (cooler colors). This visual distribution accurately reflects the extreme class imbalance and the inherent clinical difficulty of opportunistic screening for these echocardiographic targets.

C: Triangular radar plot of cross-cohort performance. The triangular radar plot compares the mean PRAUC of four models (Random Init-R18, Merl-R18, ECGCLIP-R18, and ECGCLIP-R34) across

three validation cohorts. The distance from the center to each vertex represents the mean PRAUC for that model. ECGCLIP models, particularly ECGCLIP-R34, achieve larger triangle areas, demonstrating robust cross-cohort generalization.

D: Circular bar plot of internal versus external performance. This circular bar chart compares the PRAUC values of all four models on the internal test set (darker bars) and the combined external test set (lighter bars, generated by merging the Shiyuan and BIDMC cohorts) across 39 ECHO tasks. Tasks are grouped by diagnostic category. For most tasks, ECGCLIP models maintain high and consistent PRAUC values between internal and external testing, confirming robust generalization to unseen patient populations.

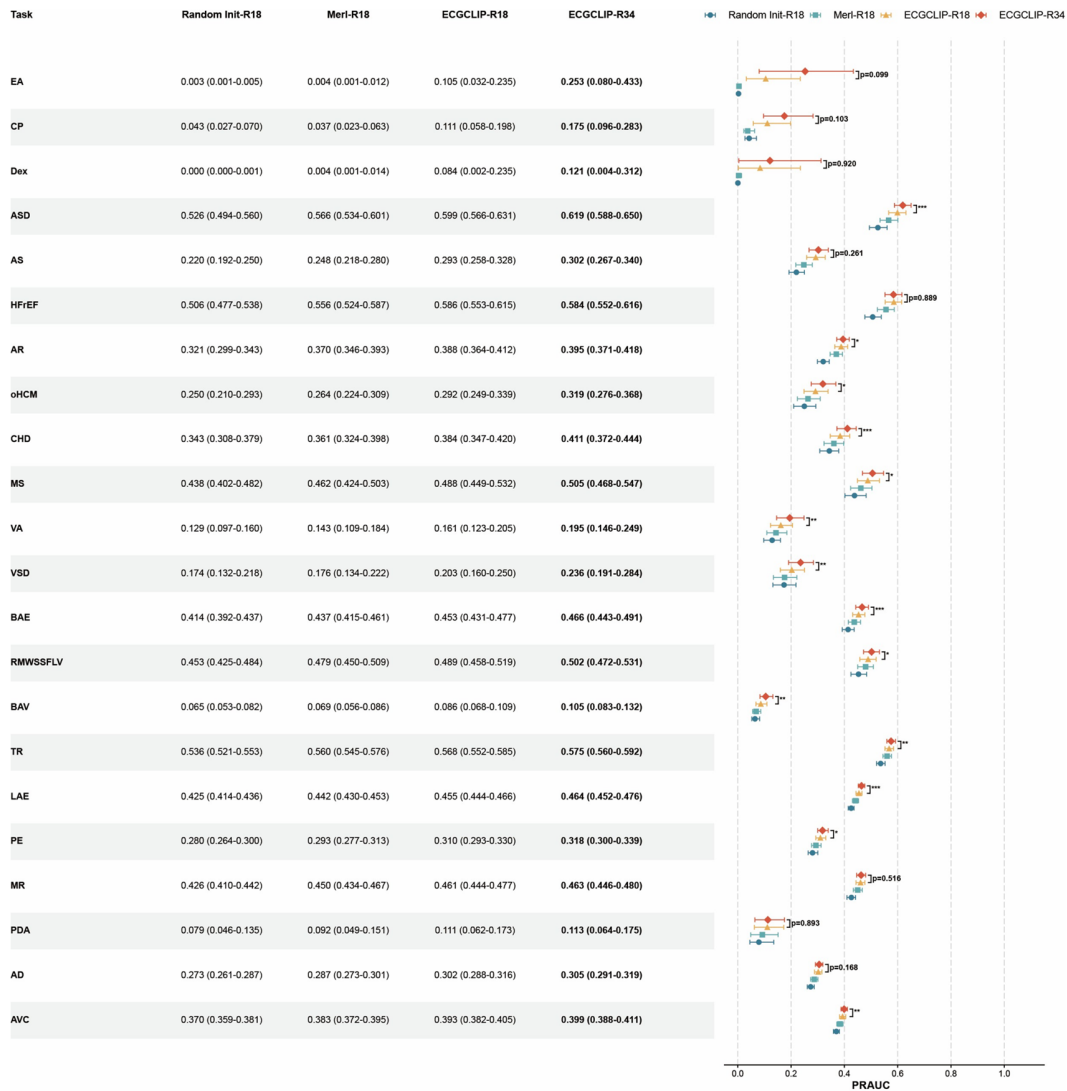


Figure 5. Performance comparison on the top 22 ECHO diagnostic tasks showing the greatest improvement.

This figure presents a combined data table (left) and forest plot (right) illustrating model performance across 22 specific echocardiographic tasks. These tasks were selected based on the greatest absolute PRAUC improvements achieved by ECGCLIP-R34 compared to the Random Init-R18 baseline. The left panel details the exact mean PRAUC alongside 95% confidence intervals (inside parentheses) for four model variants: Random Init-R18, Merl-R18, ECGCLIP-R18, and ECGCLIP-R34. The corresponding right panel visualizes these metrics, where points represent the mean PRAUC and error bars indicate the 95% confidence intervals derived from 1,000 bootstrap resamples. ECGCLIP models demonstrate substantial diagnostic improvements, particularly for low-prevalence conditions (for example, Ebstein anomaly, constrictive pericarditis, and dextrocardia), where they achieve significant PRAUC gains over near-chance baseline performance. Statistical significance between ECGCLIP-R34 and the second-best performing model for each task is indicated adjacent to the error bars by exact P values or asterisks (* $P < 0.05$, ** $P < 0.01$, *** $P < 0.001$).

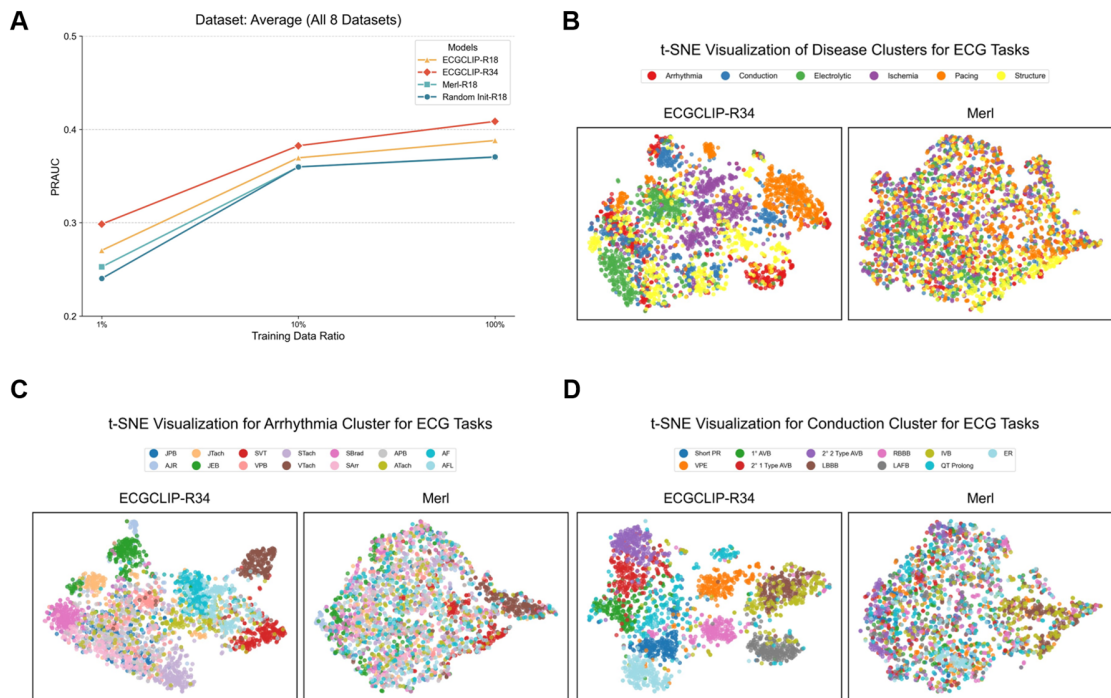


Figure 6. Data efficiency and representational quality of ECGCLIP models.

A: Data efficiency curve. This line plot compares the mean PRAUC of four models (Random Init-R18, Merl-R18, ECGCLIP-R18, and ECGCLIP-R34) across all eight ECG cohorts, evaluated at training data fractions of 1%, 10%, and 100%. ECGCLIP models, particularly ECGCLIP-R34, demonstrate high sample efficiency, matching or exceeding the performance of baseline models trained on the full dataset using only 10% of the data.

B: t-SNE visualization of global disease clusters. The t-SNE plots visualize the high-dimensional latent embeddings learned by ECGCLIP-R34 (left) and Merl-R18 (right) across 45 ECG tasks. Points are colored by six major diagnostic categories (arrhythmia, conduction, electrolytic, ischemia, pacing, and structure). ECGCLIP-R34 yields distinct, well-separated clusters, whereas Merl-R18 embeddings show heavy inter-class overlap.

C: t-SNE visualization of arrhythmia sub-clusters. This analysis focuses on arrhythmia tasks, comparing the ability of ECGCLIP-R34 (left) and Merl-R18 (right) to separate specific arrhythmia phenotypes (for example, AF, VPB, and VTach). ECGCLIP-R34 forms compact, distinct sub-clusters, while Merl-R18 representations are heavily mixed.

D: t-SNE visualization of conduction sub-clusters. This analysis compares the ability of ECGCLIP-R34 (left) and Merl-R18 (right) to separate specific conduction phenotypes (for example, RBBB, LBBB, and AVB). ECGCLIP-R34 forms clear, well-separated sub-clusters, while Merl-R18 representations exhibit heavy overlap.

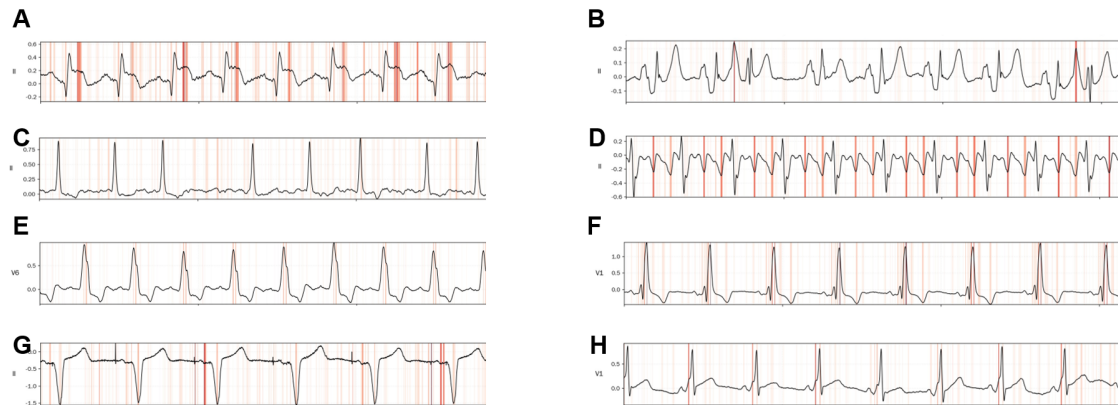


Figure 7. Interpretability of ECGCLIP predictions via Integrated Gradients (IG) saliency heatmaps.

This figure displays ECG waveforms from representative leads with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions ($IG > 0$) of the ECG signal that drive the diagnostic predictions of the model, aligning with established clinical electrocardiographic criteria.

A, ST-segment elevation myocardial infarction (STEMI): The model focuses on ST-segment elevation patterns.

B, Atrial premature beats (APC): Attention is directed to premature P-waves and associated QRS-T complexes.

C, Atrial fibrillation (AF): The model identifies global rhythm disturbances consistent with clinical criteria.

D, Atrial flutter (AFL): Characteristic sawtooth flutter waves are highlighted.

E, Left bundle branch block (LBBB): Attention aligns with the defining QRS morphology of the block.

F, Right bundle branch block (RBBB): The model focuses on diagnostic QRS patterns.

G, DDD pacing: Key highlighted regions include pacing spikes and abnormal QRS complexes.

H, Pre-excitation (WPW): The model highlights delta waves in the early QRS complex.

These saliency maps visually confirm that ECGCLIP's diagnostic decisions are based on physiologically coherent and clinically meaningful waveform features.

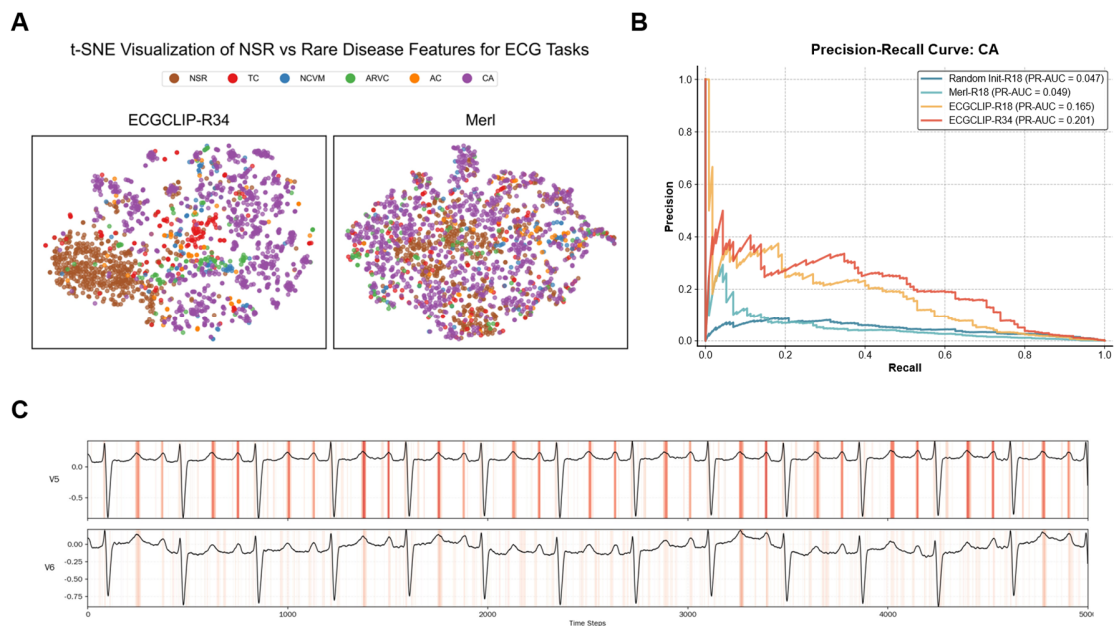


Figure 8. Diagnostic performance and interpretability of ECGCLIP for rare cardiac diseases.

A, t-SNE visualization of rare disease clusters. The t-SNE plots visualize the high-dimensional latent embeddings learned by ECGCLIP-R34 (left) and Merl-R18 (right) for distinguishing normal sinus rhythm (NSR, brown) from five rare cardiomyopathies: Takotsubo cardiomyopathy (TC, red), noncompaction of ventricular myocardium (NCVM, blue), arrhythmogenic right ventricular cardiomyopathy (ARVC, green), alcoholic cardiomyopathy (AC, orange), and cardiac amyloidosis (CA, purple). ECGCLIP-R34 yields distinct, well-separated clusters, whereas Merl-R18 embeddings show heavy inter-class overlap.

B, Precision-Recall (PR) curve for cardiac amyloidosis (CA). The PR curve compares the performance of four models (Random Init-R18, Merl-R18, ECGCLIP-R18, and ECGCLIP-R34) for the detection of CA. The PRAUC is reported in the legend. ECGCLIP models, particularly ECGCLIP-R34, achieve a substantially higher PRAUC (0.201) compared to near-random baseline performance (PRAUC < 0.05), demonstrating an enhanced capability to detect this rare and diagnostically elusive condition.

C, Integrated Gradients (IG) saliency heatmaps for CA. This figure displays ECG waveforms from leads V5 and V6 with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model prediction of CA, focusing on P wave and T wave segments that correlate with known diastolic dysfunction criteria in cardiac amyloidosis.

Supplementary Materials

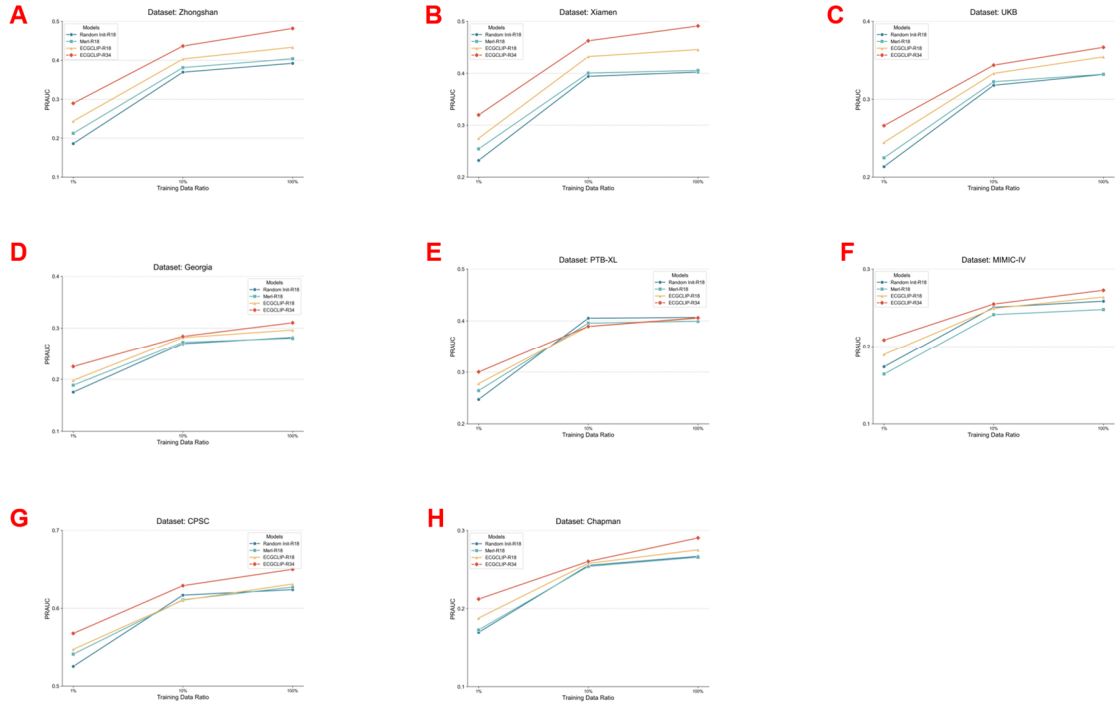


Figure S1. Data efficiency of ECGCLIP models across eight independent electrocardiographic cohorts.

This figure consists of eight line plots (A–H), each comparing the mean PRAUC of four models (Random Init-R18, Merl-R18, ECGCLIP-R18, ECGCLIP-R34) across three training data fractions (1%, 10%, and 100%) within a distinct cohort: A, Zhongshan; B, Xiamen; C, UKB; D, Georgia; E, PTB-XL; F, MIMIC-IV; G, CPSC; H, Chapman. ECGCLIP models, particularly the deeper ECGCLIP-R34 architecture, generally demonstrate superior sample efficiency. Notably, ECGCLIP-R34 matches or exceeds the performance of baseline models trained on the full dataset using only 10% of the training data in most cohorts. This superiority holds true across individual cohorts with the specific exception of the PTB-XL dataset at higher data fractions. These results highlight the exceptional data efficiency and transfer learning capability of the model across diverse patient populations and recording devices.

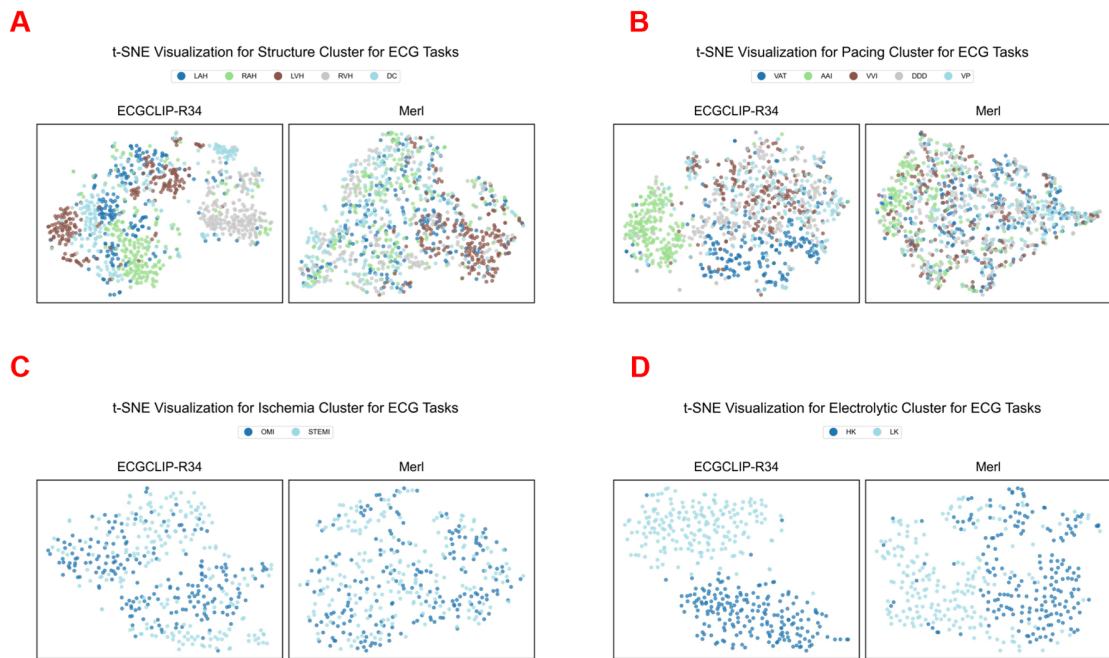


Figure S2. t-SNE visualization of sub-clusters for structure, pacing, ischemia, and electrolytic ECG tasks.

A, t-SNE visualization for the structure cluster. This analysis focuses on structure-related tasks, comparing the ability of ECGCLIP-R34 (left) and Merl (right) to separate specific structural phenotypes (e.g., LAH, RAH, LVH, RVH, DC). ECGCLIP-R34 forms compact, distinct sub-clusters, while Merl representations are heavily mixed.

B, t-SNE visualization for the pacing cluster. This analysis focuses on pacing-related tasks, comparing the ability of ECGCLIP-R34 (left) and Merl (right) to separate specific pacing phenotypes (e.g., VAT, AAI, VVI, DDD, VP). ECGCLIP-R34 forms clear, well-separated sub-clusters, while Merl representations are heavily overlapping.

C, t-SNE visualization for the ischemia cluster. This analysis focuses on ischemia-related tasks, comparing the ability of ECGCLIP-R34 (left) and Merl (right) to separate old myocardial infarction (OMI) and ST-segment elevation myocardial infarction (STEMI). ECGCLIP-R34 yields more distinct clusters, while Merl embeddings show heavy inter-class overlap.

D, t-SNE visualization for the electrolytic cluster. This analysis focuses on electrolytic tasks, comparing the ability of ECGCLIP-R34 (left) and Merl (right) to separate hyperkalemia (HK) and hypokalemia (LK). ECGCLIP-R34 forms clear, well-separated clusters, while Merl representations are heavily overlapping.

Across all sub-cluster analyses, ECGCLIP-R34 consistently learns more clinically meaningful and discriminative feature representations compared to Merl.

t-SNE Visualization of Disease Clusters for ECHO Tasks

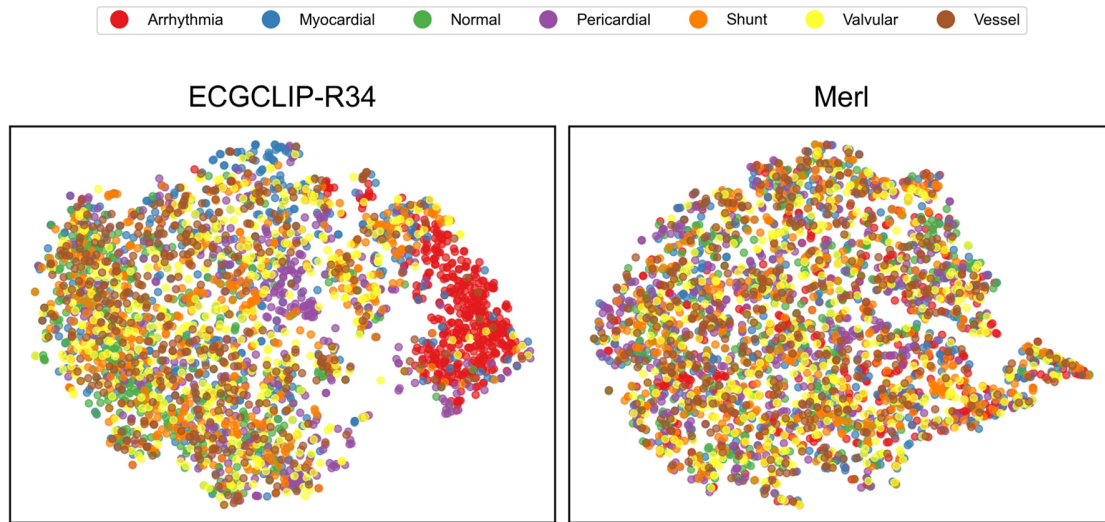


Figure S3. t-SNE visualization of disease clusters for echocardiographic (ECHO) tasks.

This figure presents t-SNE visualizations of the high-dimensional latent embeddings learned by ECGCLIP-R34 (left) and Merl (right) for 39 echocardiographic tasks. Points are colored by seven major diagnostic categories: arrhythmia (red), myocardial (blue), normal (green), pericardial (purple), shunt (orange), valvular (yellow), and vessel (brown). Given the inherent clinical difficulty of inferring structural and functional phenotypes solely from electrical signals, the clusters naturally exhibit more overlap than primary electrocardiographic tasks. Nevertheless, ECGCLIP-R34 demonstrates discernible clustering tendencies and improved regional grouping of related diagnostic categories. In contrast, MERL-R18 embeddings remain heavily entangled with no clear structural organization. This analysis confirms that while opportunistic echocardiographic screening remains challenging, ECGCLIP-R34 learns comparatively more coherent feature representations than baseline models.

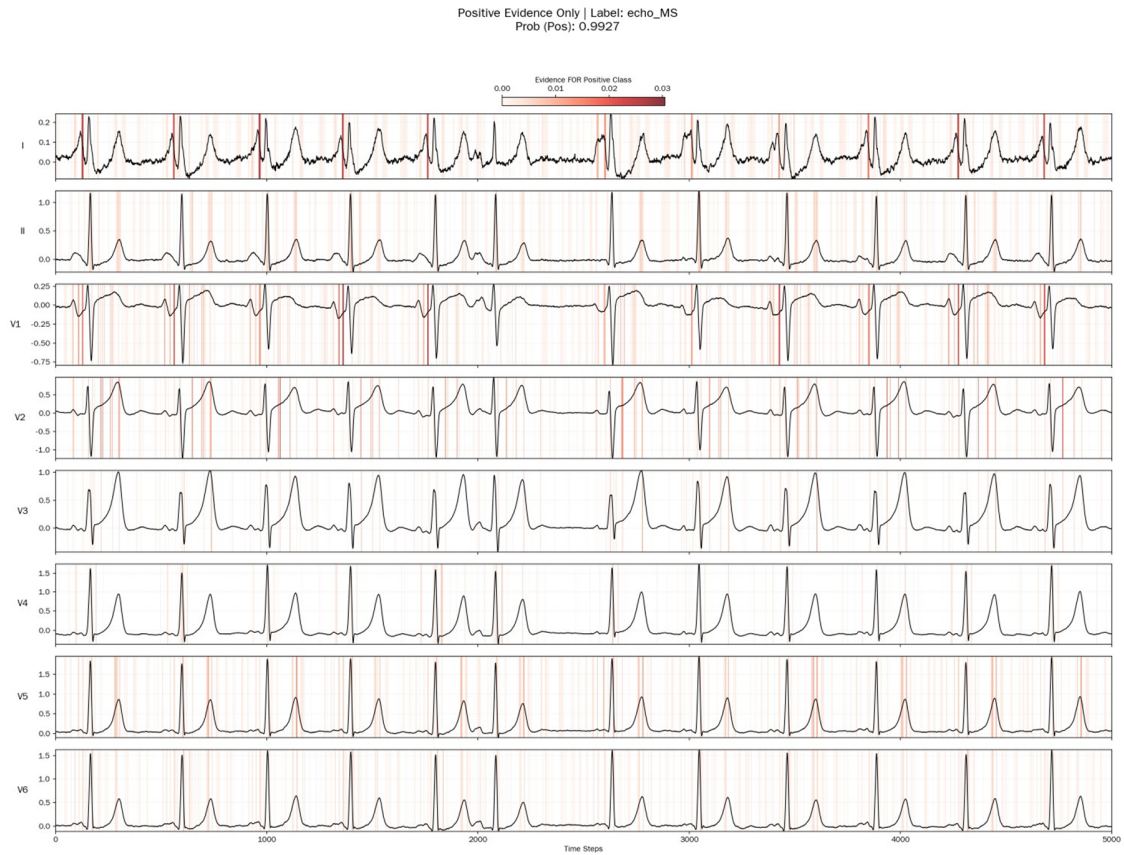


Figure S4. Integrated Gradients (IG) saliency heatmap for mitral stenosis (MS) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions ($IG > 0$) of the ECG signal that drive the model's prediction of mitral stenosis (MS), with a predicted probability of 0.9927. The model's attention is focused on waveform segments consistent with known electrocardiographic criteria for MS, including P-wave abnormalities indicative of left atrial hypertrophy and associated QRS-T complex changes.

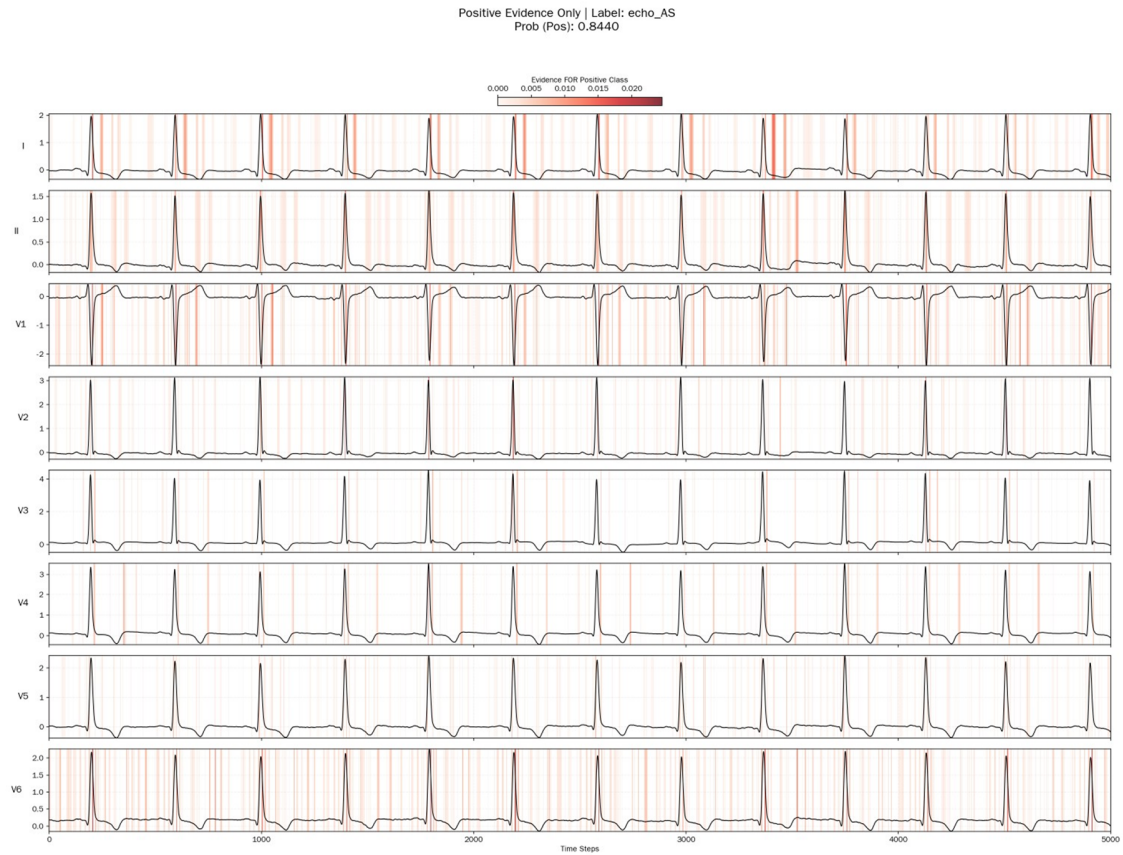


Figure S5. Integrated Gradients (IG) saliency heatmap for aortic stenosis (AS) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model's prediction of aortic stenosis (AS), with a predicted probability of 0.8440. The model's attention is focused on waveform segments consistent with known electrocardiographic criteria for AS, including left ventricular hypertrophy patterns (e.g., increased QRS amplitude, ST-T secondary changes) associated with chronic pressure overload.

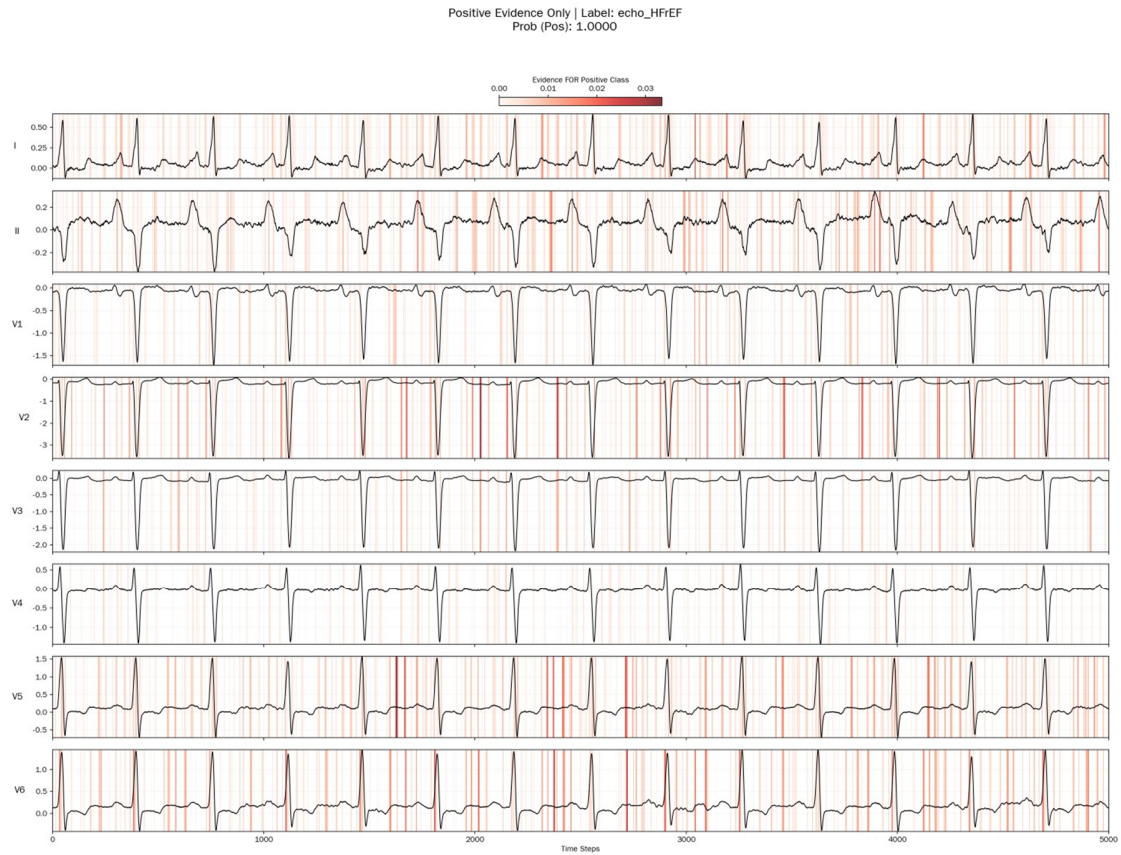


Figure S6. Integrated Gradients (IG) saliency heatmap for heart failure with reduced ejection fraction (HFrEF) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions ($IG > 0$) of the ECG signal that drive the model's prediction of heart failure with reduced ejection fraction (HFrEF), with a predicted probability of 1.0000. The model's attention is focused on waveform segments consistent with known electrocardiographic criteria for HFrEF, including low QRS voltage, ST-T abnormalities, and rhythm disturbances associated with reduced myocardial contractility and ventricular remodeling.

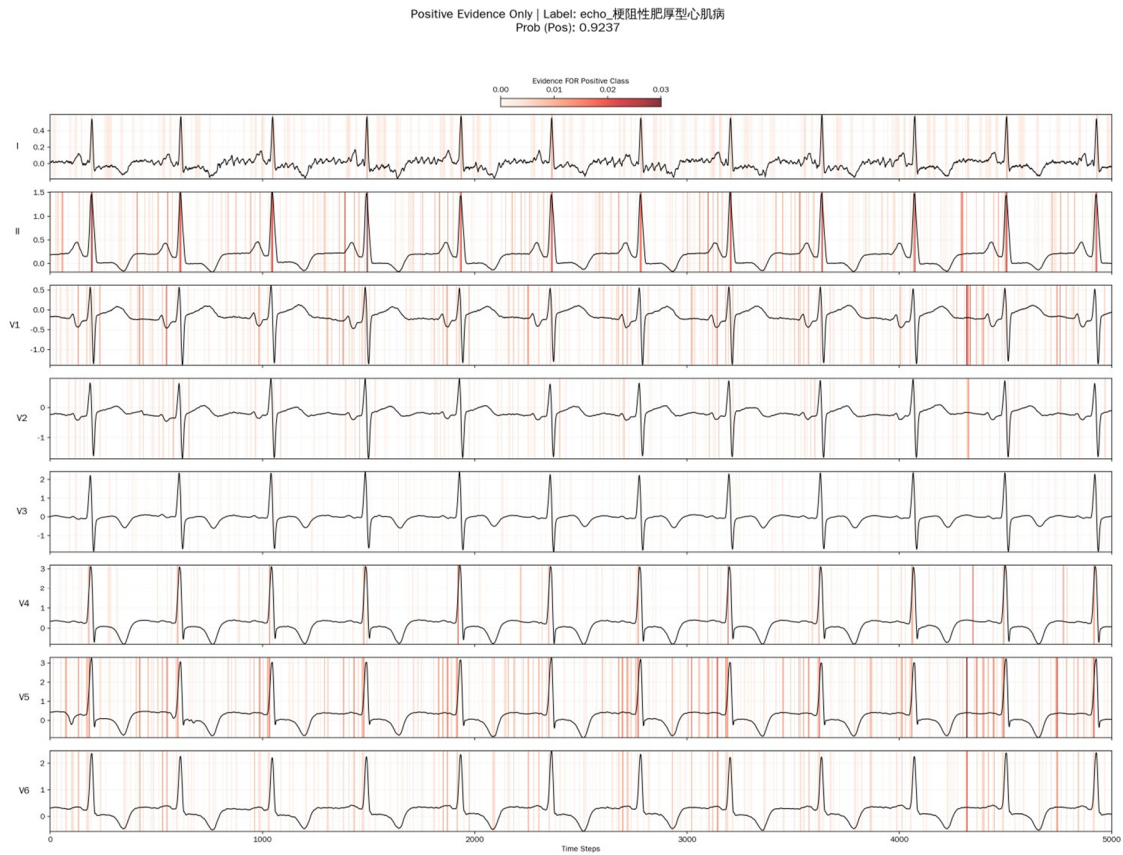


Figure S7. Integrated Gradients (IG) saliency heatmap for obstructive hypertrophic cardiomyopathy (oHCM) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1-V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model's prediction of obstructive hypertrophic cardiomyopathy (oHCM), with a predicted probability of 0.9237. The model's attention is focused on waveform segments consistent with known electrocardiographic criteria for oHCM, including left ventricular hypertrophy patterns (e.g., increased QRS amplitude, ST-T secondary changes) and associated P-wave abnormalities.

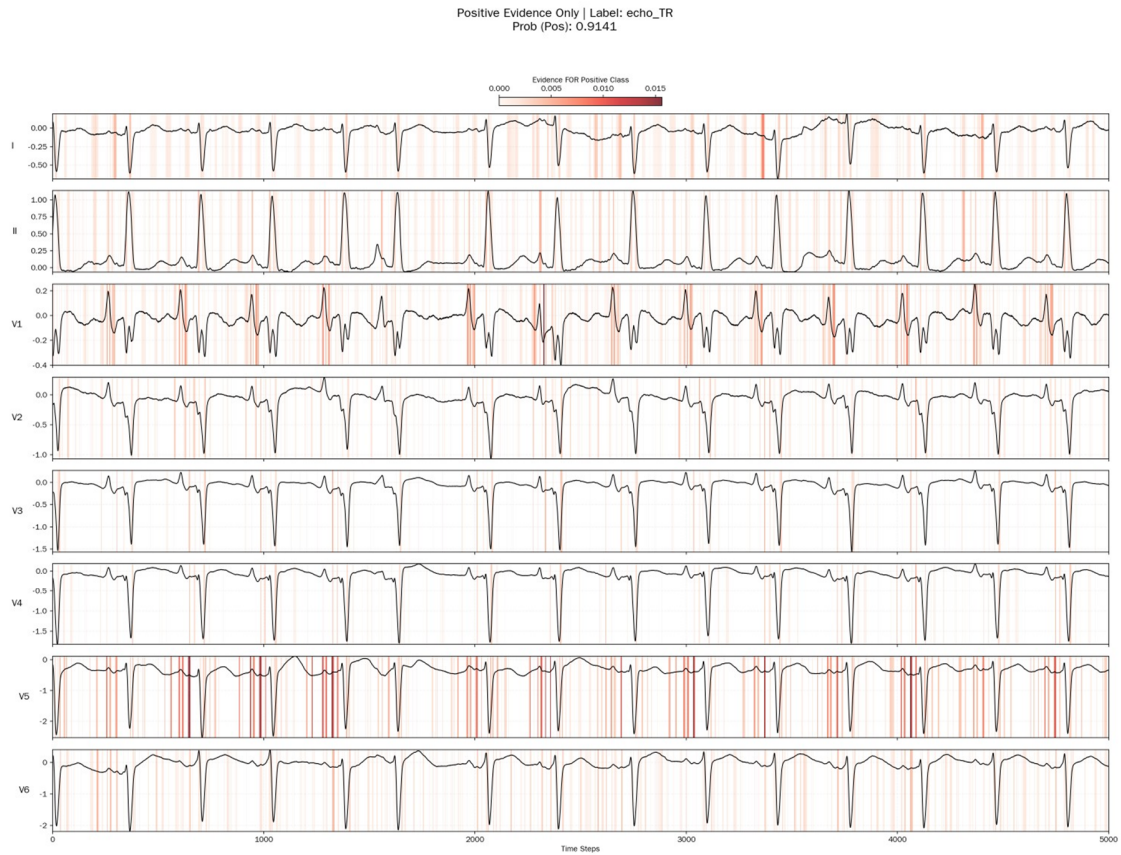


Figure S8. Integrated Gradients (IG) saliency heatmap for tricuspid regurgitation (TR) prediction. This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions ($IG > 0$) of the ECG signal that drive the model’s prediction of tricuspid regurgitation (TR), with a predicted probability of 0.9141. The model’s attention is focused on waveform segments consistent with known electrocardiographic criteria for TR, including P-wave and QRS-T complex changes indicative of right atrial and ventricular enlargement associated with chronic tricuspid regurgitation.



Figure S9. Integrated Gradients (IG) saliency heatmap for first-degree atrioventricular block (1° AVB) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model’s prediction of first-degree atrioventricular block (1° AVB), with a predicted probability of 0.9467. The model’s attention is focused on the PR segment and P-wave regions, consistent with the electrocardiographic criterion of a prolonged PR interval (> 0.20 s) for 1° AVB.

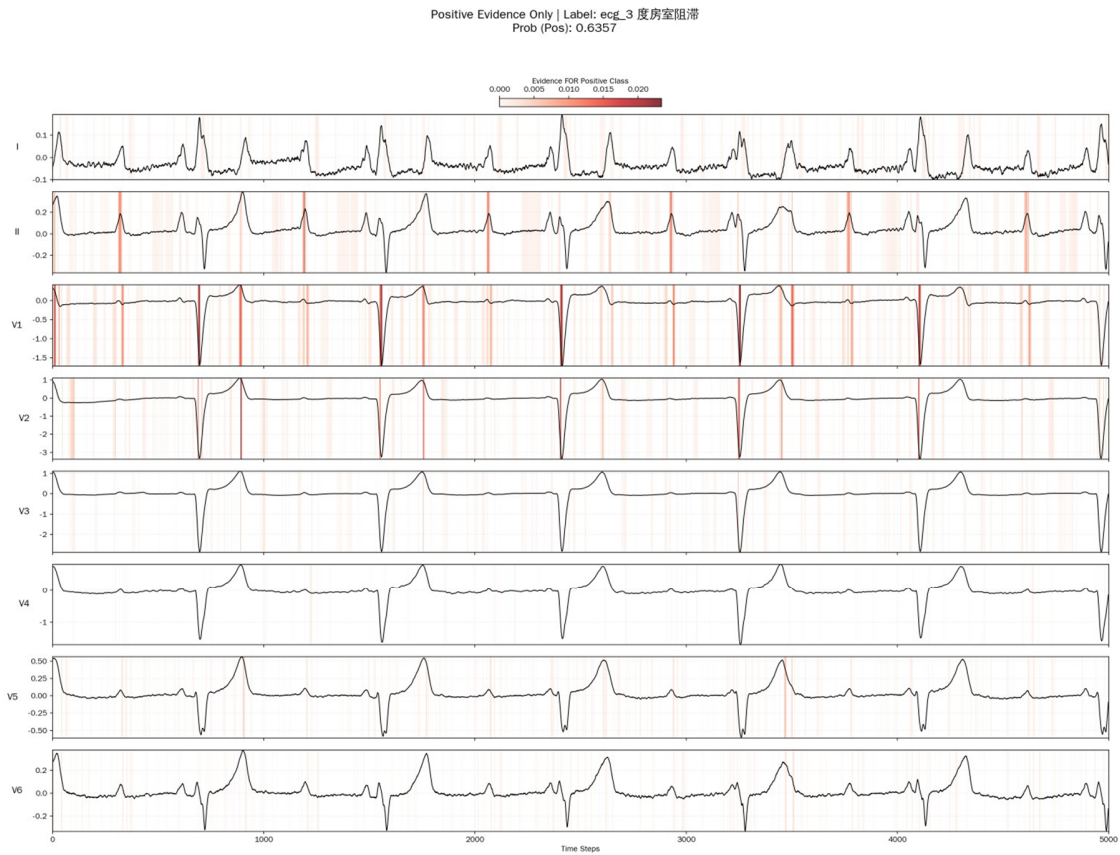


Figure S10. Integrated Gradients (IG) saliency heatmap for third-degree atrioventricular block (3° AVB) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model's prediction of third-degree atrioventricular block (3° AVB), with a predicted probability of 0.6357. The model's attention is focused on waveform segments consistent with atrioventricular (AV) dissociation, a hallmark of 3° AVB, where P waves and QRS complexes occur independently.

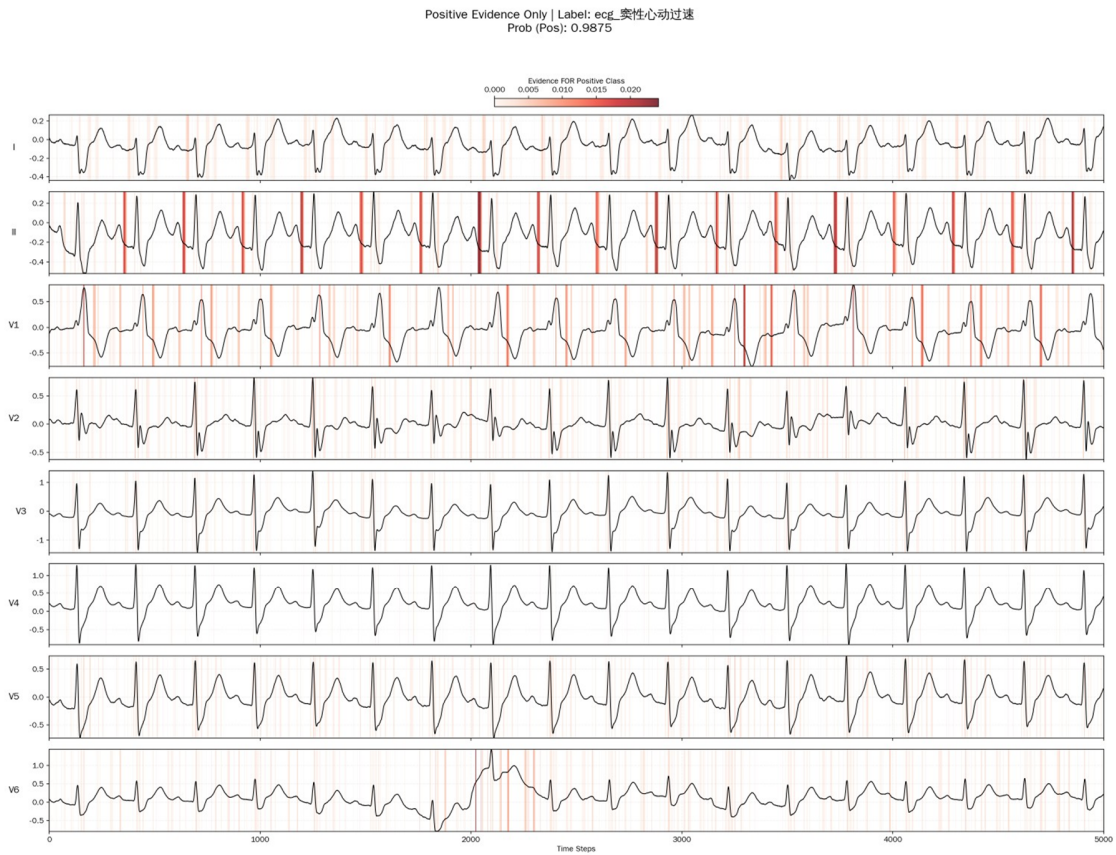


Figure S11. Integrated Gradients (IG) saliency heatmap for sinus tachycardia prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model’s prediction of sinus tachycardia, with a predicted probability of 0.9875. The model’s attention is focused on the RR intervals and overall rhythm, consistent with the electrocardiographic criterion of a heart rate > 100 bpm with sinus P waves.



Figure S12. Integrated Gradients (IG) saliency heatmap for sinus bradycardia prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model's prediction of sinus bradycardia, with a predicted probability of 0.9854. The model's attention is focused on the RR intervals and P-wave segments, consistent with the electrocardiographic criterion of a heart rate < 60 bpm with sinus P waves.

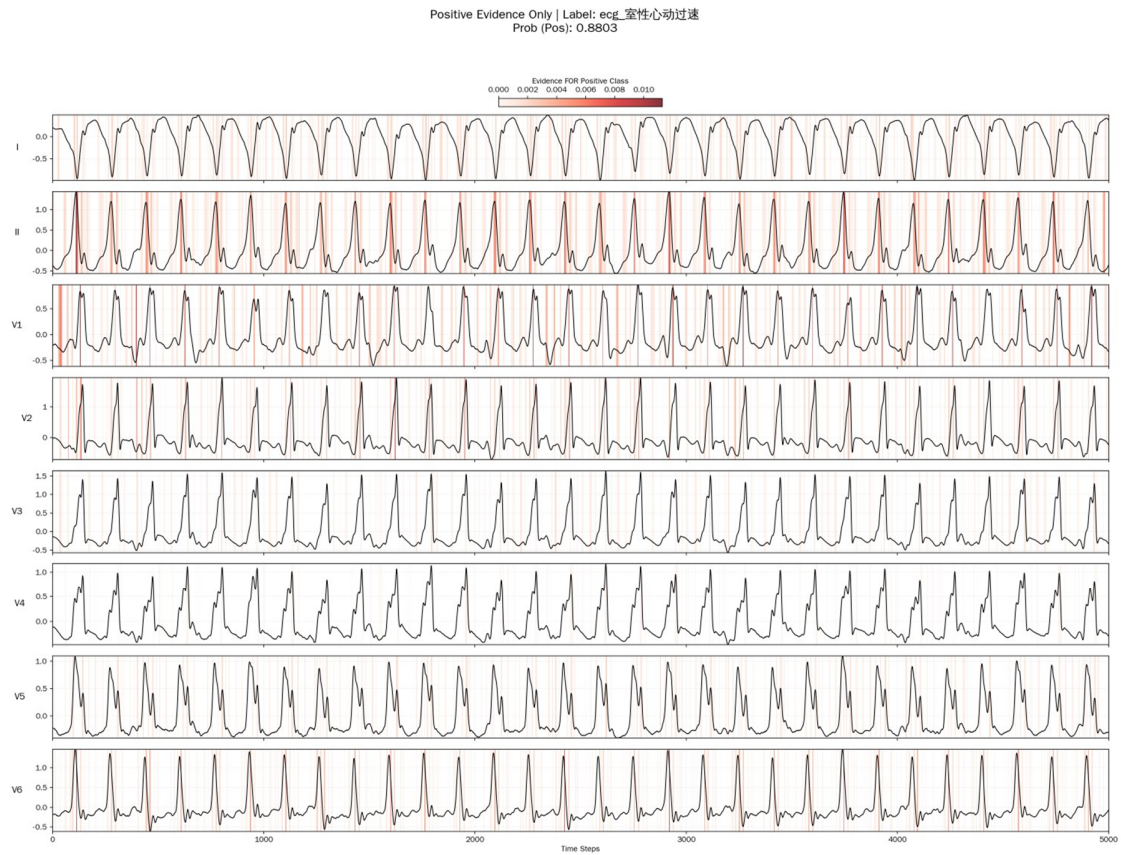


Figure S13. Integrated Gradients (IG) saliency heatmap for ventricular tachycardia (VT) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model’s prediction of ventricular tachycardia (VT), with a predicted probability of 0.8803. The model’s attention is focused on waveform segments consistent with known electrocardiographic criteria for VT, including wide, aberrant QRS complexes and rapid ventricular rhythm.

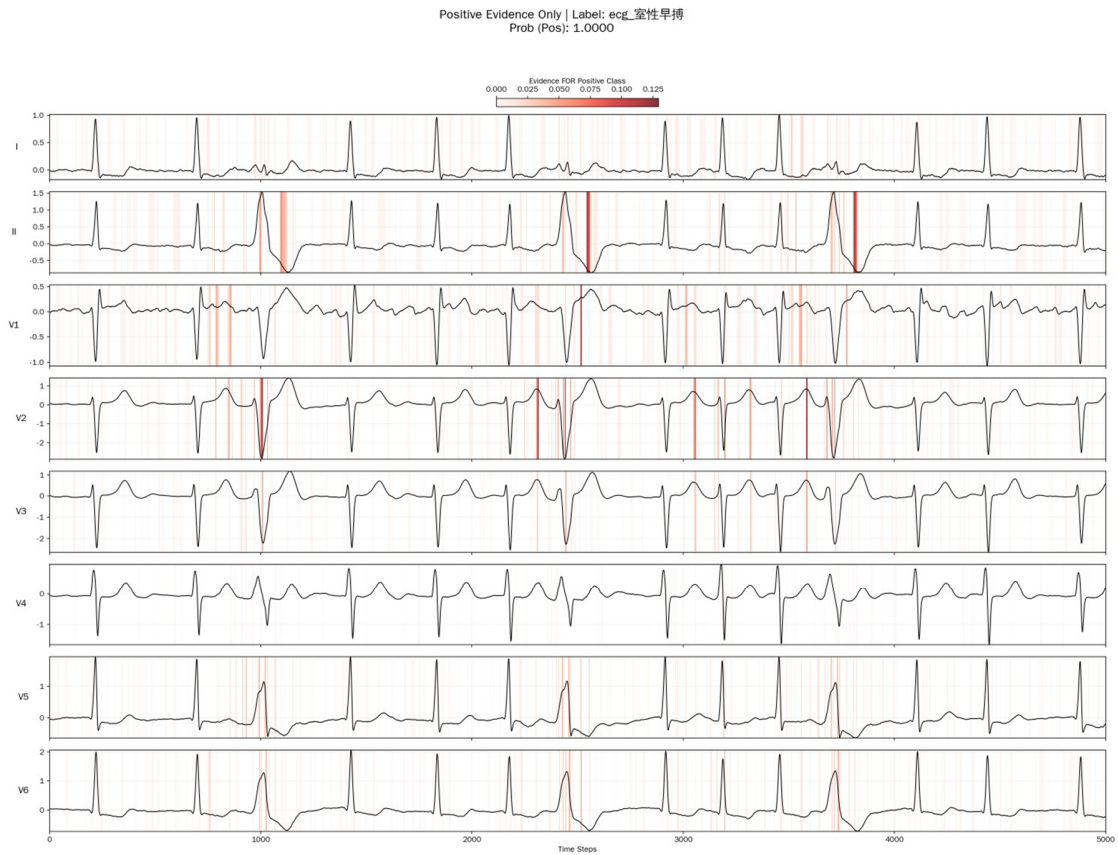


Figure S14. Integrated Gradients (IG) saliency heatmap for premature ventricular contractions (PVCs) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions ($IG > 0$) of the ECG signal that drive the model’s prediction of premature ventricular contractions (PVCs), with a predicted probability of 1.0000. The model’s attention is focused on waveform segments consistent with known electrocardiographic criteria for PVCs, including wide, premature QRS complexes with compensatory pauses.

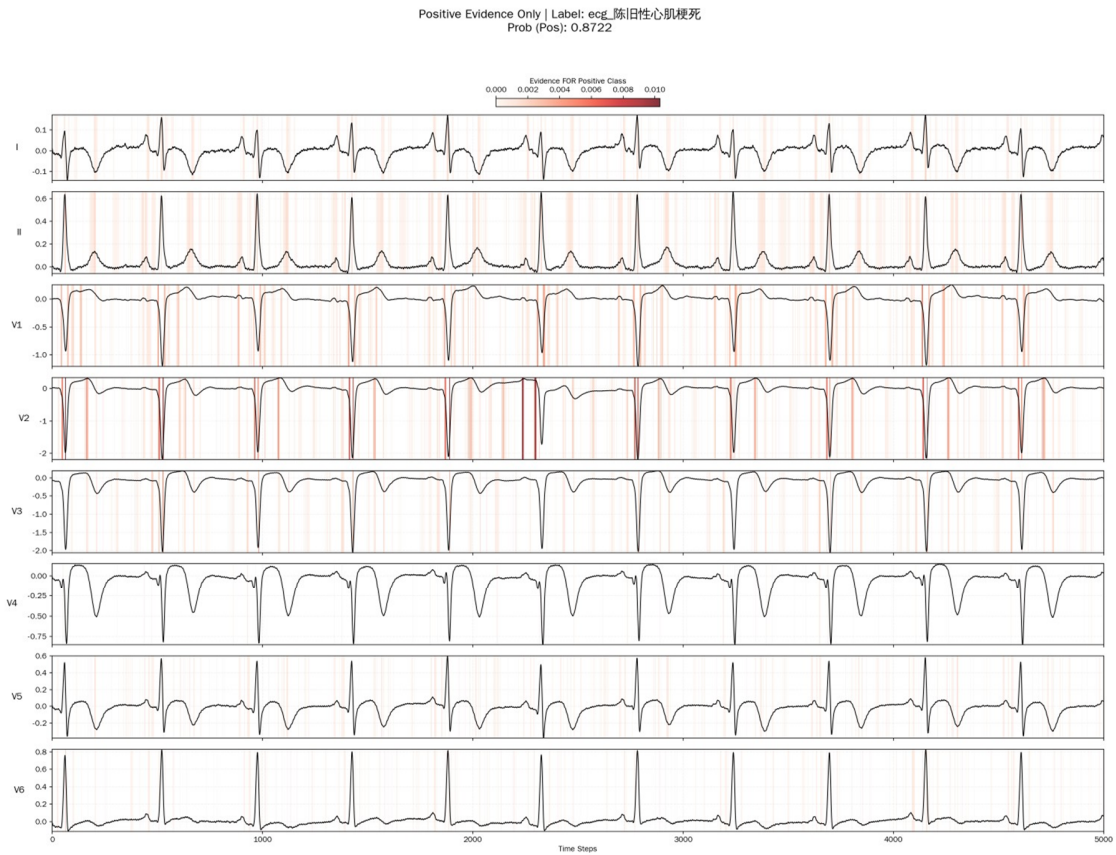


Figure S15. Integrated Gradients (IG) saliency heatmap for old myocardial infarction (OMI) prediction.

This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions (IG > 0) of the ECG signal that drive the model’s prediction of old myocardial infarction (OMI), with a predicted probability of 0.8722. The model’s attention is focused on waveform segments consistent with known electrocardiographic criteria for OMI, including pathological Q waves and associated ST-T changes indicative of prior myocardial scar.



Figure S16. Integrated Gradients (IG) saliency heatmap for atrial septal defect (ASD) prediction. This figure displays 8-lead ECG waveforms (I, II, V1–V6) with superimposed red heatmaps, generated using the Integrated Gradients method. The heatmaps highlight the positive attribution regions ($IG > 0$) of the ECG signal that drive the model’s prediction of atrial septal defect (ASD), with a predicted probability of 0.9545. The model’s attention is focused on waveform segments consistent with known electrocardiographic criteria for ASD, including P-wave and QRS-T complex changes indicative of right atrial and ventricular enlargement associated with chronic left-to-right shunting.

Table S1: Overview of data cohorts.

Cohort Category	Region (Country)	Cohort Name	# Patients	# ECG	Source
Internal Development Cohort (Pre-training & Tuning)	Asia (China)	Zhongshan	1,558,654	3,559,279	Hospital
	Asia (China)	Shiyuan	143,078	219,428	Hospital
External Test cohorts (Independent Testing)	North America (USA)	BIDMC	67,827	146,001	Hospital
	Asia (China)	Xiamen	30,032	44,362	Hospital
	North America (USA)	MIMIC-IV-ECG	160,000	772,155	Public Repository ¹
	Europe (United Kingdom)	UKB	50,407	50,407	Public Repository ²
	Asia (China)	Chapman	45,151	45,151	Public Repository ³
	Europe (Germany)	PTB-XL	21,799	21,799	Public Repository ⁴
	North America (USA)	Georgia	10,344	10,344	Public Repository ⁵
	Asia (China)	CPSC2018	6,877	6,877	Public Repository ⁶
Total	-	-	2,094,169	4,875,803	-

¹ <https://physionet.org/content/mimic-iv-ecg/1.0/>² <https://www.ukbiobank.ac.uk/>³ <https://physionet.org/content/ecg-arrhythmia/1.0.0/>⁴ <https://physionet.org/content/ptb-xl/1.0.3/>⁵ <https://physionet.org/content/challenge-2020/1.0.2/#files>⁶ <http://2018.icbeb.org/Challenge.html>

Table S2: Detailed data partitioning of the Internal Development Cohort.

Phase	Subset	Allocation	# Patients	# ECG
Pre-training	Training	98%	1,298,359	2,781,289
	Validation	2%	26,497	56,673
	Total	100%	1,324,856	2,837,962
Downstream Development	Training	70%	163,700	505,171
	Validation	15%	35,049	107,973
	Test	15%	35,049	108,173
	Total	100%	233,798	721,317

All data partitioning was strictly performed at the patient level to ensure independence between subsets and prevent data leakage.

Table S3: Definitions of diagnostic tasks and prevalence in the Internal Development Cohort.

Category	Abbreviation	Full Diagnostic Name	Zhongshan (Training)	Zhongshan (Validation)	Zhongshan (Test)
Tier1: ECG					
Interpretation					
Ischemia	OMI	Old Myocardial Infarction	5,749 (1.14%)	1,281 (1.19%)	1,411 (1.30%)
	STEMI	Acute ST-Elevation Myocardial Infarction	1,714 (0.34%)	402 (0.37%)	401 (0.37%)
Arrhythmia	NSR	Normal Sinus Rhythm	352,913 (69.86%)	75,840 (70.24%)	75,675 (69.96%)
	SBrad	Sinus Bradycardia	43,646 (8.64%)	9,408 (8.71%)	9,506 (8.79%)
	STach	Sinus Tachycardia	21,343 (4.22%)	4,592 (4.25%)	4,450 (4.11%)
	VPB	Ventricular Premature Beat	31,471 (6.23%)	6,815 (6.31%)	6,837 (6.32%)
	APB	Atrial Premature Beat	24,913 (4.93%)	5,416 (5.02%)	5,308 (4.91%)
	SArr	Sinus Arrhythmia	9,621 (1.90%)	2,134 (1.98%)	2,119 (1.96%)
	AF	Atrial Fibrillation	38,074 (7.54%)	7,769 (7.20%)	7,919 (7.32%)
	AFL	Atrial Flutter	8,328 (1.65%)	1,714 (1.59%)	1,853 (1.71%)
	JPB	Atrioventricular Junctional Premature Beat	2,667 (0.53%)	572 (0.53%)	564 (0.52%)
	AJR	Accelerated Atrioventricular Junctional Rhythm	3,454 (0.68%)	792 (0.73%)	769 (0.71%)
ATach	Atrial Tachycardia	3,074 (0.61%)	666 (0.62%)	649 (0.60%)	
JTach	Atrioventricular Junctional Tachycardia	1,521 (0.30%)	326 (0.30%)	328 (0.30%)	
JEB	Atrioventricular Junctional Escape Beat	1,727 (0.34%)	334 (0.31%)	353 (0.33%)	
SVT	Supraventricular Tachycardia	664 (0.13%)	160 (0.15%)	126 (0.12%)	
VEB	Ventricular Escape Beat	622 (0.12%)	118 (0.11%)	161 (0.15%)	
VTach	Ventricular Tachycardia	619 (0.12%)	139 (0.13%)	149 (0.14%)	
Conduction	RBBB	Right Bundle Branch Block	32,175 (6.37%)	6,991 (6.47%)	6,832 (6.32%)
Conduction	1° AVB	First-Degree Atrioventricular Block	18,671 (3.70%)	4,266 (3.95%)	3,823 (3.53%)
	QT Prolong	Prolonged QT Interval	5,963 (1.18%)	1,255 (1.16%)	1,198 (1.11%)

Category	Abbreviation	Full Diagnostic Name	Zhongshan (Training)	Zhongshan (Validation)	Zhongshan (Test)
	ER	Early Repolarization	2,328 (0.46%)	457 (0.42%)	459 (0.42%)
	LAFB	Left Anterior Fascicular Block	3,895 (0.77%)	878 (0.81%)	807 (0.75%)
	LBBB	Left Bundle Branch Block	4,481 (0.89%)	1,056 (0.98%)	828 (0.77%)
	IVB	Intraventricular Block	3,601 (0.71%)	888 (0.82%)	835 (0.77%)
	Short PR	Short PR Interval	1,010 (0.20%)	195 (0.18%)	192 (0.18%)
	VPE	Ventricular Pre-Excitation	901 (0.18%)	184 (0.17%)	202 (0.19%)
	3° AVB	Third-Degree Atrioventricular Block	1,232 (0.24%)	200 (0.19%)	249 (0.23%)
	2° 1 Type AVB	Second-Degree AV Block, Mobitz Type I	675 (0.13%)	112 (0.10%)	136 (0.13%)
	2° 2 Type AVB	Second-Degree AV Block, Mobitz Type II	299 (0.06%)	55 (0.05%)	50 (0.05%)
	LVH	Left Ventricular Hypertrophy	29,285 (5.80%)	6,528 (6.05%)	6,189 (5.72%)
	LAH	Left Atrial Hypertrophy	9,110 (1.80%)	2,058 (1.91%)	1,883 (1.74%)
	RVH	Right Ventricular Hypertrophy	3,879 (0.77%)	827 (0.77%)	984 (0.91%)
Structure	VA	Ventricular Aneurysm	888 (0.18%)	212 (0.20%)	174 (0.16%)
	RAH	Right Atrial Hypertrophy	407 (0.08%)	87 (0.08%)	89 (0.08%)
	DC	Dextrocardia	143 (0.03%)	36 (0.03%)	24 (0.02%)
	BVH	Biventricular Hypertrophy	186 (0.04%)	36 (0.03%)	44 (0.04%)
	VVI	VVI Pacing	4,338 (0.86%)	798 (0.74%)	874 (0.81%)
	VAT	VAT Pacing	2,219 (0.44%)	446 (0.41%)	443 (0.41%)
Pacing	DDD	DDD Pacing	1,013 (0.20%)	244 (0.23%)	219 (0.20%)
	AAI	AAI Pacing	875 (0.17%)	174 (0.16%)	196 (0.18%)
	VP	Ventricular Pacing	477 (0.09%)	107 (0.10%)	117 (0.11%)
Pacing	PSM	Pacing/Sensing Malfunction	193 (0.04%)	28 (0.03%)	35 (0.03%)
Electrolytic	HK	Hyperkalemia	202 (0.04%)	43 (0.04%)	45 (0.04%)

Category	Abbreviation	Full Diagnostic Name	Zhongshan (Training)	Zhongshan (Validation)	Zhongshan (Test)	
Electrolytic	LK	Hyperkalemia	213 (0.04%)	47 (0.04%)	56 (0.05%)	
Tier 2: ECHO						
Diagnosis						
Normal	NER	Normal Echocardiogram at Rest	97,681 (19.34%)	20,776 (19.24%)	21,120 (19.52%)	
	LAE	Left Atrial Enlargement	38,427 (7.61%)	8,218 (7.61%)	8,108 (7.50%)	
	HCM	Hypertrophic Cardiomyopathy	27,086 (5.36%)	5,695 (5.27%)	5,538 (5.12%)	
	nHCM	Non-Obstructive Hypertrophic Cardiomyopathy	25,126 (4.97%)	5,258 (4.87%)	5,130 (4.74%)	
	BAE	Bi-Atrial Enlargement	10,584 (2.10%)	2,258 (2.09%)	2,295 (2.12%)	
	HFrEF	Heart Failure with Reduced Ejection Fraction	5,712 (1.13%)	1,176 (1.09%)	1,242 (1.15%)	
	RMWSSFLV	Reduced Multi-Wall Segment Systolic Function of Left Ventricle	6,045 (1.20%)	1,325 (1.23%)	1,264 (1.17%)	
	Myocardial	DCM	Dilated Cardiomyopathy	3,216 (0.64%)	646 (0.60%)	622 (0.58%)
		oHCM	Obstructive Hypertrophic Cardiomyopathy	2,005 (0.40%)	452 (0.42%)	418 (0.39%)
		ICM	Ischemic Cardiomyopathy	1,363 (0.27%)	334 (0.31%)	319 (0.29%)
		LVDD	Left Ventricular Diastolic Dysfunction	1,929 (0.38%)	431 (0.40%)	445 (0.41%)
		VA	Ventricular Aneurysm	606 (0.12%)	164 (0.15%)	175 (0.16%)
RAE		Right Atrial Enlargement	444 (0.09%)	88 (0.08%)	105 (0.10%)	
AM		Atrial Myxoma	322 (0.06%)	60 (0.06%)	81 (0.07%)	
LAM		Left Atrial Mass	155 (0.03%)	44 (0.04%)	37 (0.03%)	
Valvular	AVC	Aortic Valve Calcification	34,312 (6.79%)	7,410 (6.86%)	7,218 (6.67%)	
	TR	Tricuspid Regurgitation	19,123 (3.79%)	3,941 (3.65%)	4,071 (3.76%)	
	MR	Mitral Regurgitation	16,375 (3.24%)	3,521 (3.26%)	3,426 (3.17%)	
	AR	Aortic Regurgitation	9,135 (1.81%)	1,935 (1.79%)	1,914 (1.77%)	
	MVPLAC	Mitral Valve Posterior Leaflet Annular	5,135 (1.02%)	1,106 (1.02%)	1,141 (1.05%)	

Category	Abbreviation	Full Diagnostic Name	Zhongshan (Training)	Zhongshan (Validation)	Zhongshan (Test)
		Calcification			
	MS	Mitral Stenosis	4,018 (0.80%)	799 (0.74%)	757 (0.70%)
	AS	Aortic Stenosis	2,983 (0.59%)	639 (0.59%)	653 (0.60%)
	BAV	Bicuspid Aortic Valve	2,304 (0.46%)	517 (0.48%)	479 (0.44%)
	PS	Pulmonary Stenosis	161 (0.03%)	40 (0.04%)	49 (0.05%)
	EA	Ebstein's Anomaly	120 (0.02%)	36 (0.03%)	27 (0.02%)
	PR	Pulmonary Regurgitation	99 (0.02%)	18 (0.02%)	21 (0.02%)
	TS	Tricuspid Stenosis	36 (0.01%)	2 (0.00%)	12 (0.01%)
	CHD	Congenital Heart Disease	4,139 (0.82%)	920 (0.85%)	838 (0.77%)
	ASD	Atrial Septal Defect	4,557 (0.90%)	1,003 (0.93%)	952 (0.88%)
	VSD	Ventricular Septal Defect	1,664 (0.33%)	359 (0.33%)	356 (0.33%)
Shunt	PFO	Patent Foramen Ovale	1,037 (0.21%)	199 (0.18%)	200 (0.18%)
	PDA	Patent Ductus Arteriosus	608 (0.12%)	119 (0.11%)	123 (0.11%)
	TOF	Tetralogy of Fallot	6 (0.00%)	2 (0.00%)	2 (0.00%)
	PH	Pulmonary Hypertension	32,557 (6.44%)	7,028 (6.51%)	6,917 (6.39%)
Vessel	AD	Aortic Dilatation	18,831 (3.73%)	4,076 (3.78%)	3,936 (3.64%)
	Dex	Dextrocardia	99 (0.02%)	15 (0.01%)	15 (0.01%)
Pericardial	PE	Pericardial Effusion	11,926 (2.36%)	2,530 (2.34%)	2,445 (2.26%)
Pericardial	CP	Constrictive Pericarditis	248 (0.05%)	68 (0.06%)	58 (0.05%)
Arrhythmia	AF	Atrial Fibrillation	64,960 (12.86%)	13,255 (12.28%)	13,615 (12.59%)
Tier 3: Rare Diseases					
	CA	Cardiac Amyloidosis	194 (0.04%)	85 (0.08%)	115 (0.11%)
-	ARVC	Arrhythmogenic Right Ventricular Cardiomyopathy	51 (0.01%)	7 (0.01%)	4 (0.00%)

Category	Abbreviation	Full Diagnostic Name	Zhongshan (Training)	Zhongshan (Validation)	Zhongshan (Test)
	TC	Takotsubo Cardiomyopathy	64 (0.01%)	17 (0.02%)	14 (0.01%)
	AC	Alcoholic Cardiomyopathy	28 (0.01%)	2 (0.00%)	6 (0.01%)
	NCVM	Noncompaction of Ventricular Myocardium	39 (0.01%)	1 (0.00%)	16 (0.01%)

Values are presented as No. (%), representing the count of positive cases and the prevalence rate within each specific data split (Training, Validation, or Test).

Table S4: Diagnostic label mapping criteria for the external Shiyuan test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria (English Translation [Original Chinese])
NER	“Resting...No abnormalities seen” (静息...未见异常)
LAE	“Left Atrium...Enlarged” (左房...大)
HCM	“Thickening” (增厚) OR “Hypertrophy” (肥厚)
nHCM	“Non-obstructive” (非梗阻)
BAE	“Bi-Atrial...Enlarged” (双房...大)
HFrEF	†
RMWSSFLV	“LV...Systolic...Weakened” (左室...收缩...减弱) AND NOT “Non-global” (非整体)
DCM	“Dilated...Heart” (扩张...心)
oHCM	“Obstructive” (梗阻) AND NOT Negation terms (e.g., "No/Not/None" [未/不/无])
LVDD	“LV Diastolic” (左室舒张)
VA	“Ventricular Aneurysm” (室壁瘤)
RAE	“Right Atrium...Enlarged” (右房...大)
AM	“Myxoma” (粘液瘤)
LAM	“Left Atrium...Occupying Lesion” (左房...占位)
AVC	“Aorta...Calcification” (主动脉...钙化)
TR	†
MR	†
AR	†
MVPLAC	“Mitral Valve...Calcification” (二尖瓣...钙化)
MS	†
AS	†
BAV	“Two-leaf...Aorta” (二叶...主动脉)
PS	†
EA	“Tricuspid downward displacement” (三尖瓣下移畸形)
PR	†

Target Label (Zhongshan Ontology)	Mapping Criteria (English Translation [Original Chinese])
CHD	“Congenital Heart Disease” (先天性心脏病)
ASD	“Atrial Septal...Defect” (房间隔...缺损)
VSD	“Ventricular Septal...Defect” (室间隔...缺损)
PFO	“Patent Foramen Ovale” (卵圆孔未闭)
PDA	“Patent Ductus Arteriosus” (动脉导管未闭)
TOF	“Fallot” (法洛)
PH	“Pulmonary Artery...High Pressure” (肺动脉...高压)
AD	“Aorta...Widened” (主动脉...增宽)
Dex	“Dextrocardia” (右位心)
PE	“Pericardial Effusion” (心包积液)
CP	“Constrictive...Pericarditis” (缩窄...心包炎)

†: Indicates labels retrieved directly from structured database fields (pre-existing diagnoses) within the Shiyuan dataset. All other labels were extracted via regular expression matching applied to the original clinical reports. Matching keywords are presented as English translations followed by the representative original Chinese terms (in parentheses).

Table S5: Diagnostic label mapping criteria for the external BIDMC test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria (Structured Variables & Thresholds)
LAE	la_size > 40 mm
BAE	ra_size AND la_size > 40 mm
HFrEF	biplane_lvef < 40 %
oHCM	lv_obstruction / lv_obstruction_valsalva
VA	lv_aneurysm
RAE	ra_size > 40 mm
TR	tricuspid_regurg
MR	mitral_regurg
AR	aortic_regurg
MS	mitral_stenosis
AS	aortic_stenosis
PS	pulm_stenosis
PR	pulm_regurg
VSD	vsd
PH	pasp_underestimated_tr > 40 mmHg
AD	ascending_diam > 37 mm
PE	pericardial_effusion

The mapping criteria involve either applying clinical thresholds to quantitative measurements (e.g., chamber sizes) or identifying positive diagnostic flags in categorical fields.

Table S6: Diagnostic label mapping criteria for the external MIMIC-IV-ECG test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
OMI	old inferior myocardial infarction / old anterolateral myocardial infarction / old anterior myocardial infarction / old septal myocardial infarction / old lateral myocardial infarction / old posterior myocardial infarction / old anteroseptal myocardial infarction
STEMI	inferior myocardial infarction / anterior myocardial infarction / anteroseptal myocardial infarction / septal myocardial infarction / acute myocardial infarction / acute ST segment elevation myocardial infarction / lateral myocardial infarction / extensive myocardial infarction / inferior ST segment elevation / lateral ST segment elevation / anterolateral myocardial infarction / anteroseptal ST segment elevation / anterolateral ST segment elevation / inferior and lateral ST segment elevation / anterior ST segment elevation / extensive ST segment elevation / extensive ST elevation / septal ST segment elevation / acute inferior myocardial infarction / acute anterior myocardial infarction / acute lateral myocardial infarction / acute anterolateral myocardial infarction / inferior and anterior/septal ST segment elevation / acute posterior myocardial infarction / acute anteroseptal myocardial infarction / inferior and septal ST segment elevation
NSR	normal ECG
SBrad	sinus bradycardia / marked sinus bradycardia / irregular sinus bradycardia / slow sinus arrhythmia
STach	sinus tachycardia / irregular sinus tachycardia
VPB	ventricular premature complex / frequent ventricular premature complex / frequent multifocal ventricular premature complex / multifocal ventricular premature complex / interpolated ventricular premature complex / bigeminal ventricular premature complex / multiple premature ventricular complexes / multiform ventricular premature complexes / paired ventricular premature complexes / ventricular bigeminy / ventricular couplets / ventricular trigeminy / ventricular pre-excitation / multifocal interpolated ventricular premature complex
APB	atrial premature complex / bigeminal atrial premature complex / frequent atrial premature complex / occasional frequent atrial premature complex / atrial premature complexes in couplets
SArr	sinus arrhythmia
AF	atrial fibrillation

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
JPB	supraventricular premature beats / supraventricular bigeminy / aberrantly conducted supraventricular premature beats
AJR	accelerated junctional rhythm
ATach	ectopic atrial tachycardia / atrial tachycardia / sinus/ectopic atrial tachycardia
JTach	junctional tachycardia
SVT	supraventricular tachycardia
VEB	paroxysmal idioventricular rhythm
VTach	nonsustained ventricular tachycardia / ventricular tachycardia / unsustained ventricular tachycardia
RBBB	right bundle branch block / incomplete right bundle branch block / atypical right bundle branch block
1° AVB	borderline first degree atrioventricular block / first degree atrioventricular block
QT Prolong	prolonged QT interval / long QTc interval
ER	early repolarization
LAFB	left anterior fascicular block / anterior fascicular block
LBBB	left bundle branch block / incomplete left bundle branch block / atypical left bundle branch block
IVB	intraventricular conduction defect / non-specific intraventricular conduction delay / non-specific intraventricular conduction disturbance (block) / extensive intraventricular conduction defect
Short PR	short PR interval
VPE	Wolf-Parkinson-White syndrome
3° AVB	complete atrioventricular block / third degree atrioventricular block
2° 1 Type AVB	second degree atrioventricular block
2° 2 Type AVB	4:1 atrioventricular block / 3:1 atrioventricular block / second degree atrioventricular block / predominant 4:1 atrioventricular block / predominant 3:1 atrioventricular block
LVH	left ventricular hypertrophy / voltage criteria for left ventricular hypertrophy
LAH	left atrial enlargement / bi-atrial enlargement
RVH	right ventricular hypertrophy / bi-atrial enlargement

Target Label	Mapping Criteria (Original Diagnostic Terms)
(Zhongshan Ontology)	
RAH	right atrial enlargement
DC	dextrocardia
BVH	ventricular hypertrophy / biventricular hypertrophy
VAT	dual chamber pacemaker
DDD	dual chamber pacemaker / atrioventricular sequential pacing pattern
HK	hyperkalemia

The mapping schema aggregates diverse original diagnostic terms and synonym variations into the target standardized categories.

Table S7: Diagnostic label mapping criteria for the external UKB test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria
STEMI	“ACUTE MI” OR “acute...infarction”
NSR	“Normal ECG” OR “Normal Sinus Rhythm”
SBrad	“sinus bradycardia”
VPB	“premature ventricular”
APB	“premature atrial complexes”
SArr	“sinus arrhythmia”
AF	“Atrial fibrillation”
AFL	“Atrial flutter”
JPB	“premature supraventricular” OR “junctional pacemaker”
AJR	“Junctional rhythm”
JTach	“Junctional rhythm”
JEB	“Junctional rhythm”
RBBB	“Right bundle branch block”
1° AVB	“1st degree AV block”
QT Prolong	“Prolonged QT”
ER	“Early repolarization”
LAFB	“Left anterior fascicular block”
LBBB	“Left bundle branch block”
IVB	“Nonspecific intraventricular”
Short PR	“short PR”
LVH	“Left ventricular hypertrophy”
LAH	“Left atrial enlargement” OR “Biatrial enlargement”
RVH	“Right ventricular hypertrophy”
RAH	Right atrial enlargement” OR “Biatrial enlargement”

All labels were extracted via regular expression matching applied to the original clinical reports.

Table S8: Diagnostic label mapping criteria for the external Chapman test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
STEMI	MI (myocardial infarction) / MIBW (myocardial infarction in back wall) / MIFW (Myocardial infarction in the front wall) / MILW (Myocardial infarction in the lower wall) / MISW (Myocardial infarction in the side wall)
SBrad	SB (Sinus Bradycardia)
STach	ST (Sinus Tachycardia)
VPB	VB (ventricular bigeminy) / VPB (ventricular premature beat) / VPE (ventricular preexcitation)
APB	ABI (atrial bigeminy) / APB (atrial premature beats)
SArr	SA (Sinus Irregularity) / SAAWR (Sinus Atrium to Atrial Wandering Rhythm)
AF	AFIB (Atrial Fibrillation)
AFL	AF (Atrial Flutter)
JPB	JPT (junctional premature beat)
ATach	AT (Atrial Tachycardia)
JEB	JEB (junctional escape beat)
SVT	SVT (Supraventricular Tachycardia) / AVNRT (Atrioventricular Node Reentrant Tachycardia) / AVRT (Atrioventricular Reentrant Tachycardia)
VEB	VEB (ventricular escape beat) / VET (ventricular escape trigeminy)
RBBB	RBBB (right bundle branch block)
1° AVB	1AVB (1 degree atrioventricular block) / AVB (atrioventricular block)
QT Prolong	QTIE (QT interval extension)
ER	ERV (Early repolarization of the ventricles)
LAFB	LFBBB (left front bundle branch block)
IVB	IVB (Intraventricular block)
VPE	WPW (Wolf-Parkinson-White syndrome)
3° AVB	3AVB (3 degree atrioventricular block) / AVB (atrioventricular block)

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
2° 1 Type AVB	2AVB (2 degree atrioventricular block) / 2AVB1 (2 degree atrioventricular block(Type one)) / AVB (atrioventricular block)
2° 2 Type AVB	2AVB (2 degree atrioventricular block) / 2AVB2 (2 degree atrioventricular block(Type two)) / AVB (atrioventricular block)
LVH	LVH (left ventricle hypertrophy)
RVH	RVH (right ventricle hypertrophy)
RAH	RAH (right atrial hypertrophy)
LK	UW (U wave)

The mapping schema aggregates diverse original diagnostic terms and synonym variations into the target standardized categories.

Table S9: Diagnostic label mapping criteria for the external PTBXL test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
STEMI	IMI (inferior myocardial infarction) / ASMI (anteroseptal myocardial infarction) / ILMI (inferolateral myocardial infarction) / AMI (anterior myocardial infarction) / ALMI (anterolateral myocardial infarction) / LMI (lateral myocardial infarction) / IPLMI (inferoposterolateral myocardial infarction) / IPMI (inferoposterior myocardial infarction) / PMI (posterior myocardial infarction)
NSR	NORM (normal ECG)
SBrad	SBRAD (sinus bradycardia)
STach	STACH (sinus tachycardia)
VPB	PVC (ventricular premature complex)
APB	PAC (atrial premature complex)
SArr	SARRH (sinus arrhythmia)
AF	AFIB (atrial fibrillation)
AFL	AFLT (atrial flutter)
AJR	SVTAC (supraventricular tachycardia)
ATach	SVTAC (supraventricular tachycardia)
JTach	SVTAC (supraventricular tachycardia)
SVT	SVARR (supraventricular arrhythmia) / PSVT (paroxysmal supraventricular tachycardia)
RBBB	IRBBB (incomplete right bundle branch block) / CRBBB (complete right bundle branch block)
1° AVB	1AVB (first degree AV block)
QT Prolong	LNGQT (long QT-interval)
LAFB	LAFB (left anterior fascicular block)
LBBS	ILBBB (incomplete left bundle branch block) / CLBBB (complete left bundle branch block)
IVB	IVCD (non-specific intraventricular conduction disturbance (block))
VPE	WPW (Wolf-Parkinson-White syndrome)
3° AVB	3AVB (third degree AV block)

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
2° 1 Type AVB	2AVB (second degree AV block)
2° 2 Type AVB	2AVB (second degree AV block)
LVH	LVH (left ventricular hypertrophy) / VCLVH (voltage criteria (QRS) for left ventricular hypertrophy)
LAH	LAO/LAE (left atrial overload/enlargement)
RVH	RVH (right ventricular hypertrophy)
VA	ANEUR (ST-T changes compatible with ventricular aneurysm)
RAH	RAO/RAE (right atrial overload/enlargement)
VVI	PACE (normal functioning artificial pacemaker)
VAT	PACE (normal functioning artificial pacemaker)
DDD	PACE (normal functioning artificial pacemaker)
AAI	PACE (normal functioning artificial pacemaker)
VP	PACE (normal functioning artificial pacemaker)

The mapping schema aggregates diverse original diagnostic terms and synonym variations into the target standardized categories.

Table S10: Diagnostic label mapping criteria for the external Georgia test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
STEMI	acute myocardial infarction / anterior myocardial infarction / myocardial infarction
SBrad	Bradycardia (Brady) / sinus bradycardia (SB)
STach	sinus tachycardia (STach)
VPB	paired ventricular premature complexes (VPVC) / premature ventricular contractions (PVC) / ventricular premature beats (VPB) / ventricular bigeminy (VBIG) / ventricular ectopics (VEB) / ventricular pre-excitation (VPEX) / ventricular trigeminy (VTrig)
APB	atrial bigeminy (AB) / premature atrial contraction (PAC)
SArr	sinus arrhythmia (SA)
AF	atrial fibrillation (AF) / chronic atrial fibrillation (CAF) / paroxysmal atrial fibrillation (PAF) / rapid atrial fibrillation (RAF)
AFL	atrial flutter (AFL)
JPB	junctional premature complex (JPC) / supraventricular bigeminy (SVB) / supraventricular premature beats (SVPB)
AJR	accelerated junctional rhythm (AJR)
ATach	atrial tachycardia (ATach)
JTach	junctional tachycardia (JTach)
JEB	junctional escape (JE)
SVT	paroxysmal supraventricular tachycardia (PSVT) / supraventricular tachycardia (SVT)
VEB	ventricular escape beat (VEsB) / ventricular escape rhythm (VEsR)
RBBB	bundle branch block (BBB) / complete right bundle branch block (CRBBB) / incomplete right bundle branch block (IRBBB) / right bundle branch block (RBBB)
1° AVB	1st degree av block (IAVB) / av block (AVB)
QT Prolong	prolonged qt interval (LQT)
ER	early repolarization (ERe)

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
LAFB	left anterior fascicular block (LAnFB)
LBBB	bundle branch block (BBB) / incomplete left bundle branch block (ILBBB) / left bundle branch block (LBBB)
IVB	diffuse intraventricular block (DIB) / nonspecific intraventricular conduction disorder (NSIVCB)
Short PR	shortened pr interval (SPRI)
VPE	wolff parkinson white pattern (WPW)
3° AVB	av block (AVB) / complete heart block (CHB)
2° 1 Type AVB	2nd degree av block (IIAVB) / av block (AVB) / mobitz type i wenckebach atrioventricular block (MoI)
2° 2 Type AVB	2nd degree av block (IIAVB) / av block (AVB)
LVH	left ventricular hypertrophy (LVH) / ventricular hypertrophy (VH)
LAH	atrial hypertrophy (AH) / left atrial enlargement (LAE) / left atrial hypertrophy (LAH)
RVH	right ventricular hypertrophy (RVH) / ventricular hypertrophy (VH)
RAH	atrial hypertrophy (AH) / right atrial hypertrophy (RAH) /
AAI	atrial pacing pattern (AP)
VP	ventricular pacing pattern (VPP)

The mapping schema aggregates diverse original diagnostic terms and synonym variations into the target standardized categories.

Table S11: Diagnostic label mapping criteria for the external CPSC2018 test cohort.

Target Label (Zhongshan Ontology)	Mapping Criteria (Original Diagnostic Terms)
STEMI	ST-segment elevated (STE)
NSR	Normal
VPB	Premature ventricular contraction (PVC)
APB	Premature atrial contraction (PAC)
AF	Atrial fibrillation (AF)
RBBB	Right bundle branch block (RBBB)
1° AVB	First-degree atrioventricular block (I-AVB)
LBBB	Left bundle branch block (LBBB)

The mapping schema aggregates diverse original diagnostic terms and synonym variations into the target standardized categories.

Table S12: Availability of diagnostic tasks across cohorts.

Cohort	Tier 1: ECG Interpretation	Tier 2: ECHO Diagnosis	Tier 3: Rare Disease Diagnosis
Zhongshan	√	√	√
Shiyuan	-	√	-
BIDMC	-	√	-
Xiamen	√	-	-
MIMIC-IV-ECG	√	-	-
UKB	√	-	-
Chapman	√	-	-
PTB-XL	√	-	-
Georgia	√	-	-
CPSC2018	√	-	-

√: Indicates that validated ground truth labels for this task tier were available and mapped to the standardized ontology. -: Indicates that the specific ground truth labels were unavailable and the cohort was not utilized for this specific tier.

Table S13: Prevalence of ECG diagnostic labels in External Test Cohorts.

Disease	Xiamen	MIMIC-IV-ECG	UKB	Chapman	PTBXL	Georgia	CPSC2018
OMI	47 (0.11%)	14,581 (1.91%)	-	-	-	-	-
STEMI	432 (0.97%)	166,317 (21.80%)	220 (0.44%)	123 (0.27%)	5,288 (24.26%)	7 (0.07%)	220 (3.20%)
NSR	34,336 (77.40%)	117,266 (15.37%)	47,890 (95.01%)	-	9,514 (43.64%)	-	918 (13.35%)
SBrad	5,396 (12.16%)	97,656 (12.80%)	22,546 (44.73%)	16,525 (36.68%)	637 (2.92%)	1,682 (16.27%)	-
STach	1,394 (3.14%)	67,450 (8.84%)	-	7,236 (16.06%)	826 (3.79%)	1,259 (12.18%)	-
VPB	1,277 (2.88%)	50,962 (6.68%)	1,409 (2.80%)	309 (0.69%)	1,143 (5.24%)	399 (3.86%)	700 (10.18%)
APB	1,003 (2.26%)	34,983 (4.59%)	1,049 (2.08%)	3 (0.01%)	398 (1.83%)	638 (6.17%)	616 (8.96%)
SArr	1,522 (3.43%)	25,668 (3.36%)	2,377 (4.72%)	2,535 (5.63%)	772 (3.54%)	454 (4.39%)	-
AF	427 (0.96%)	78,363 (10.27%)	772 (1.53%)	1,780 (3.95%)	1,514 (6.95%)	570 (5.51%)	1,221 (17.75%)
AFL	114 (0.26%)	-	65 (0.13%)	8,046 (17.86%)	73 (0.33%)	186 (1.80%)	-
JPB	250 (0.56%)	2,195 (0.29%)	453 (0.90%)	11 (0.02%)	-	1 (0.01%)	-
AJR	-	3,030 (0.40%)	184 (0.37%)	-	27 (0.12%)	19 (0.18%)	-
ATach	145 (0.33%)	1,884 (0.25%)	-	296 (0.66%)	27 (0.12%)	28 (0.27%)	-
JTach	114 (0.26%)	168 (0.02%)	184 (0.37%)	-	27 (0.12%)	4 (0.04%)	-
JEB	16 (0.04%)	-	184 (0.37%)	75 (0.17%)	-	5 (0.05%)	-
SVT	16 (0.04%)	939 (0.12%)	-	749 (1.66%)	179 (0.82%)	32 (0.31%)	-
VEB	15 (0.03%)	583 (0.08%)	-	63 (0.14%)	-	1 (0.01%)	-
VTach	11 (0.02%)	675 (0.09%)	-	-	-	-	-
RBBB	1,274 (2.87%)	63,641 (8.34%)	2,519 (5.00%)	649 (1.44%)	1,658 (7.61%)	992 (9.60%)	1,857 (27.00%)
1° AVB	985 (2.22%)	57,497 (7.54%)	2,954 (5.86%)	1,380 (3.06%)	793 (3.64%)	843 (8.15%)	722 (10.50%)
QT Prolong	780 (1.76%)	38,296 (5.02%)	510 (1.01%)	394 (0.87%)	117 (0.54%)	1,390 (13.45%)	-
ER	366 (0.83%)	3,912 (0.51%)	435 (0.86%)	365 (0.81%)	-	140 (1.35%)	-
LAFB	98 (0.22%)	47,563 (6.23%)	756 (1.50%)	240 (0.53%)	1,623 (7.45%)	180 (1.74%)	-
LBFB	139 (0.31%)	28,535 (3.74%)	501 (0.99%)	-	613 (2.81%)	327 (3.16%)	236 (3.43%)
IVB	57 (0.13%)	31,683 (4.15%)	308 (0.61%)	770 (1.71%)	787 (3.61%)	203 (1.96%)	-

Disease	Xiamen	MIMIC-IV-ECG	UKB	Chapman	PTBXL	Georgia	CPSC2018
Short PR	80 (0.18%)	11,012 (1.44%)	437 (0.87%)	-	-	2 (0.02%)	-
VPE	82 (0.18%)	499 (0.07%)	-	72 (0.16%)	79 (0.36%)	2 (0.02%)	-
3° AVB	2 (0.00%)	204 (0.03%)	-	320 (0.71%)	16 (0.07%)	82 (0.79%)	-
2° 1 Type AVB	11 (0.02%)	298 (0.04%)	-	339 (0.75%)	14 (0.06%)	96 (0.93%)	-
2° 2 Type AVB	5 (0.01%)	1,534 (0.20%)	-	308 (0.68%)	14 (0.06%)	96 (0.93%)	-
LVH	1,781 (4.01%)	71,267 (9.34%)	1,033 (2.05%)	645 (1.43%)	2,354 (10.80%)	1,264 (12.23%)	-
LAH	1,225 (2.76%)	16,785 (2.20%)	3,278 (6.50%)	-	426 (1.95%)	930 (9.00%)	-
RVH	192 (0.43%)	7,910 (1.04%)	88 (0.17%)	110 (0.24%)	126 (0.58%)	71 (0.69%)	-
VA	19 (0.04%)	-	-	-	104 (0.48%)	-	-
RAH	49 (0.11%)	330 (0.04%)	112 (0.22%)	36 (0.08%)	99 (0.45%)	141 (1.36%)	-
DC	7 (0.02%)	218 (0.03%)	-	-	-	-	-
BVH	-	17,772 (2.33%)	-	-	-	-	-
VVI	44 (0.10%)	-	-	-	294 (1.35%)	-	-
VAT	32 (0.07%)	239 (0.03%)	-	-	294 (1.35%)	-	-
DDD	3 (0.01%)	288 (0.04%)	-	-	294 (1.35%)	-	-
AAI	25 (0.06%)	-	-	-	294 (1.35%)	52 (0.50%)	-
VP	-	-	-	-	294 (1.35%)	45 (0.44%)	-
PSM	3 (0.01%)	-	-	-	-	-	-
HK	5 (0.01%)	1,650 (0.22%)	-	-	-	-	-
LK	-	-	-	136 (0.30%)	-	-	-

Data are expressed as No. (%), denoting the number of identified positive cases and their prevalence within the specific cohort. -: Indicates the label was not available in the source dataset.

Table S14: Prevalence of ECHO diagnostic labels in External Test Cohorts.

Disease	Shiyuan	BIDMC
NER	440 (0.20%)	-
LAE	25,834 (11.78%)	69,199 (52.43%)
HCM	17,740 (8.09%)	-
nHCM	17,258 (7.87%)	-
BAE	3,631 (1.66%)	39,303 (35.27%)
HFrEF	4,005 (1.83%)	2,385 (10.31%)
RMWSSFLV	5 (0.00%)	-
DCM	3 (0.00%)	-
oHCM	565 (0.26%)	3,208 (2.81%)
ICM	-	-
LVDD	148,641 (67.80%)	-
VA	756 (0.34%)	682 (89.50%)
RAE	218 (0.10%)	50,221 (44.46%)
AM	30 (0.01%)	-
LAM	55 (0.03%)	-
AVC	23,654 (10.79%)	-
TR	7,382 (3.40%)	15,537 (11.49%)
MR	5,244 (2.42%)	15,062 (11.16%)
AR	3,539 (1.63%)	3,000 (5.52%)
MVPLAC	3,322 (1.52%)	-
MS	287 (0.13%)	500 (2.14%)
AS	547 (0.25%)	7,266 (43.10%)
BAV	249 (0.11%)	-
PS	31 (0.01%)	70 (0.18%)
EA	10 (0.00%)	-

Disease	Shiyuan	BIDMC
PR	163 (0.08%)	2,040 (1.99%)
TS	-	-
CHD	746 (0.34%)	-
ASD	657 (0.30%)	-
VSD	139 (0.06%)	254 (73.62%)
PFO	697 (0.32%)	-
PDA	65 (0.03%)	-
TOF	27 (0.01%)	-
PH	26,621 (12.14%)	15,988 (14.61%)
AD	27,990 (12.77%)	19,699 (16.66%)
Dex	78 (0.04%)	-
PE	11,795 (5.38%)	1,677 (5.85%)
CP	4 (0.00%)	-
AF	-	-

Data are expressed as No. (%), denoting the number of identified positive cases and their prevalence within the specific cohort. -: Indicates the label was not available in the source dataset.

Table S15: Global performance overview (PRAUC) of comparative models across diverse diagnostic tiers and cohorts.

Cohort	Random Init-	Random Init-	Random Init-	ECGCLIP-	ECGCLIP-	ECGCLIP-	Merl-R18
	R18	R34	R50	R18	R34	R50	
Tier1: ECG							
Interpretation							
Zhongshan	0.392±0.325	0.427±0.322	0.404±0.313	0.433±0.324	0.482±0.326	<u>0.473±0.316</u>	0.404±0.326
Xiamen	0.402±0.349	0.444±0.339	0.415±0.337	0.446±0.340	<u>0.491±0.338</u>	0.496±0.329	0.405±0.346
MIMIC-IV-ECG	0.259±0.309	<u>0.267±0.311</u>	0.256±0.306	0.264±0.310	0.273±0.315	0.264±0.307	0.248±0.299
UKB	0.332±0.354	0.345±0.366	0.328±0.352	<u>0.354±0.362</u>	0.367±0.366	0.349±0.361	0.332±0.356
Chapman	0.267±0.302	0.274±0.306	0.266±0.305	0.275±0.303	0.290±0.304	<u>0.287±0.305</u>	0.266±0.301
PTBXL	<u>0.407±0.323</u>	0.402±0.321	0.396±0.310	0.405±0.326	0.406±0.324	0.407±0.319	0.399±0.334
Georgia	0.282±0.307	0.294±0.313	0.284±0.307	0.296±0.310	0.310±0.315	<u>0.304±0.312</u>	0.280±0.311
CPSC2018	0.624±0.335	0.637±0.330	0.619±0.337	0.631±0.330	0.650±0.313	<u>0.641±0.323</u>	0.627±0.334
Average	0.370	0.386	0.371	0.388	0.409	<u>0.403</u>	0.370
Tier 2: ECHO							
Diagnosis							
Zhongshan	0.263±0.237	0.262±0.234	0.263±0.220	0.292±0.242	<u>0.306±0.241</u>	0.310±0.241	0.274±0.245
Shiyuan	0.163±0.194	0.161±0.193	0.163±0.188	0.178±0.202	<u>0.185±0.204</u>	0.186±0.201	0.170±0.198
BIDMC	0.369±0.283	0.364±0.287	0.365±0.280	0.378±0.285	<u>0.384±0.288</u>	0.384±0.283	0.375±0.284
Average	0.265	0.263	0.264	0.283	<u>0.292</u>	0.293	0.273
Tier 3: Rare							
Disease Diagnosis							
Zhongshan	0.016	0.011	0.022	<u>0.035</u>	0.045	0.028	0.011

Reported values represent the mean ± standard deviation, averaged across all available diagnostic tasks within each specific dataset and tier. **Bold** denotes the highest PRAUC score, and underlines denotes the second-highest score within each dataset row.

Table S16: Training hyperparameters for ECGCLIP-R34 pre-training.

Hyperparameter	Value
Global Batch Size	1,400(700 per GPU × 2 GPUs)
Optimizer	AdamW
AdamW β	(0.9,0.999)
Weight Decay	1e-5
Peak Learning Rate	2e-4
Learning Rate Schedule	Cosine Annealing with Restarts ($T_0=5k$ steps, $\eta_{min}=1e-8$)
Epochs	20
Precision Type	bfloat 16 (AMP)
CLIP Temperature	0.07
Dropout (UMA module)	0.1
Text Encoder Frozen	First 9 layers (of 12)

Experiments were conducted on two NVIDIA A100 80GB GPUs, utilizing a local batch size of 700 per GPU with mixed-precision training for enhanced efficiency. Training proceeded for up to 20 epochs, with the best-performing model checkpoint selected based on the highest cross-modal retrieval recall observed on the pre-training dataset’s validation set.

Table S17: Training hyperparameters for downstream fine-tuning.

Hyperparameter	Value
Batchsize	256
Optimizer	AdamW
AdamW β	(0.9,0.999)
Weight Decay	5e-2
Peak Learning Rate	5e-4
Learning Rate Schedule	Cosine
Iterations	4,000

Performance was monitored every 200 iterations, and the checkpoint achieving the highest PRAUC on the validation set was selected for final evaluation on the independent test sets.

Table S18: Detailed task-specific performance comparison on the internal Zhongshan test cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
Tier1: ECG							
Interpretation							
OMI	Random Init-R18	0.425 (0.396 – 0.452)	0.976 (0.973 – 0.978)	0.508 (0.482 – 0.535)	0.992 (0.992 – 0.993)	0.482 (0.460 – 0.503)	0.476 (0.453 – 0.497)
	ECGCLIP-R18	0.483 (0.454 – 0.512)	0.981 (0.978 – 0.983)	0.601 (0.576 – 0.628)	0.990 (0.989 – 0.991)	0.511 (0.490 – 0.531)	0.509 (0.489 – 0.530)
	ECGCLIP-R34	0.506 (0.478 – 0.535)	0.983 (0.980 – 0.985)	0.597 (0.571 – 0.623)	0.991 (0.990 – 0.991)	0.521 (0.500 – 0.543)	0.519 (0.497 – 0.540)
	Merl-R18	0.458 (0.431 – 0.485)	0.980 (0.977 – 0.982)	0.554 (0.530 – 0.581)	0.991 (0.991 – 0.992)	0.499 (0.479 – 0.521)	0.494 (0.474 – 0.516)
STEMI	Random Init-R18	0.199 (0.164 – 0.240)	0.960 (0.951 – 0.969)	0.324 (0.279 – 0.369)	0.996 (0.996 – 0.997)	0.282 (0.243 – 0.320)	0.281 (0.242 – 0.320)
	ECGCLIP-R18	0.265 (0.223 – 0.313)	0.976 (0.969 – 0.982)	0.419 (0.370 – 0.468)	0.996 (0.996 – 0.997)	0.344 (0.306 – 0.384)	0.347 (0.308 – 0.387)
	ECGCLIP-R34	0.383 (0.332 – 0.436)	0.984 (0.978 – 0.988)	0.596 (0.550 – 0.643)	0.995 (0.995 – 0.996)	0.419 (0.384 – 0.457)	0.436 (0.401 – 0.473)
	Merl-R18	0.206 (0.169 – 0.248)	0.965 (0.957 – 0.972)	0.416 (0.366 – 0.464)	0.994 (0.994 – 0.995)	0.281 (0.247 – 0.313)	0.293 (0.258 – 0.327)
NSR	Random Init-R18	0.985 (0.984 – 0.986)	0.974 (0.973 – 0.976)	0.972 (0.971 – 0.973)	0.876 (0.873 – 0.880)	0.960 (0.959 – 0.961)	0.863 (0.860 – 0.866)
	ECGCLIP-R18	0.987 (0.986 – 0.988)	0.978 (0.977 – 0.979)	0.971 (0.969 – 0.972)	0.905 (0.901 – 0.908)	0.965 (0.964 – 0.966)	0.882 (0.879 – 0.886)
	ECGCLIP-R34	0.989 (0.988 – 0.990)	0.982 (0.981 – 0.983)	0.981 (0.981 – 0.982)	0.909 (0.905 – 0.911)	0.971 (0.970 – 0.972)	0.903 (0.900 – 0.906)
	Merl-R18	0.985 (0.985 – 0.986)	0.975 (0.974 – 0.976)	0.967 (0.966 – 0.969)	0.892 (0.888 – 0.895)	0.961 (0.960 – 0.962)	0.868 (0.864 – 0.871)
SBrad	Random Init-R18	0.927 (0.921 – 0.932)	0.991 (0.990 – 0.992)	0.877 (0.870 – 0.883)	0.986 (0.985 – 0.986)	0.866 (0.860 – 0.871)	0.853 (0.847 – 0.858)
	ECGCLIP-R18	0.942 (0.937 – 0.946)	0.993 (0.992 – 0.993)	0.883 (0.877 – 0.890)	0.988 (0.988 – 0.989)	0.882 (0.876 – 0.886)	0.870 (0.864 – 0.875)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STach	ECGCLIP-R34	0.961 (0.957 - 0.964)	0.995 (0.994 - 0.996)	0.914 (0.909 - 0.919)	0.991 (0.991 - 0.992)	0.912 (0.908 - 0.917)	0.904 (0.899 - 0.908)
	Merl-R18	0.927 (0.922 - 0.932)	0.991 (0.990 - 0.992)	0.869 (0.862 - 0.876)	0.986 (0.986 - 0.987)	0.864 (0.859 - 0.869)	0.851 (0.845 - 0.857)
	Random Init-R18	0.943 (0.937 - 0.948)	0.995 (0.994 - 0.996)	0.864 (0.854 - 0.874)	0.996 (0.996 - 0.996)	0.884 (0.877 - 0.891)	0.879 (0.872 - 0.886)
	ECGCLIP-R18	0.953 (0.948 - 0.958)	0.996 (0.996 - 0.997)	0.882 (0.871 - 0.892)	0.996 (0.996 - 0.997)	0.896 (0.889 - 0.903)	0.892 (0.885 - 0.899)
	ECGCLIP-R34	0.967 (0.962 - 0.971)	0.997 (0.997 - 0.998)	0.920 (0.912 - 0.928)	0.997 (0.996 - 0.997)	0.922 (0.917 - 0.928)	0.919 (0.913 - 0.925)
	Merl-R18	0.943 (0.938 - 0.949)	0.995 (0.994 - 0.996)	0.875 (0.865 - 0.885)	0.995 (0.995 - 0.996)	0.883 (0.876 - 0.890)	0.878 (0.871 - 0.885)
VPB	Random Init-R18	0.768 (0.757 - 0.778)	0.956 (0.954 - 0.959)	0.671 (0.660 - 0.682)	0.989 (0.988 - 0.990)	0.733 (0.724 - 0.742)	0.720 (0.711 - 0.730)
	ECGCLIP-R18	0.833 (0.824 - 0.841)	0.970 (0.967 - 0.972)	0.730 (0.719 - 0.740)	0.992 (0.991 - 0.993)	0.790 (0.782 - 0.798)	0.780 (0.772 - 0.788)
	ECGCLIP-R34	0.913 (0.907 - 0.919)	0.983 (0.981 - 0.985)	0.846 (0.837 - 0.854)	0.994 (0.994 - 0.995)	0.875 (0.869 - 0.881)	0.867 (0.861 - 0.874)
	Merl-R18	0.795 (0.786 - 0.804)	0.962 (0.959 - 0.964)	0.682 (0.671 - 0.693)	0.991 (0.991 - 0.992)	0.753 (0.744 - 0.761)	0.742 (0.733 - 0.751)
APB	Random Init-R18	0.255 (0.244 - 0.266)	0.862 (0.857 - 0.866)	0.486 (0.472 - 0.498)	0.928 (0.926 - 0.929)	0.337 (0.327 - 0.346)	0.308 (0.298 - 0.318)
	ECGCLIP-R18	0.343 (0.330 - 0.355)	0.889 (0.885 - 0.893)	0.464 (0.451 - 0.477)	0.956 (0.955 - 0.957)	0.401 (0.390 - 0.412)	0.369 (0.358 - 0.380)
	ECGCLIP-R34	0.589 (0.574 - 0.604)	0.937 (0.933 - 0.940)	0.563 (0.549 - 0.576)	0.982 (0.981 - 0.982)	0.587 (0.575 - 0.598)	0.567 (0.555 - 0.579)
SArr	Merl-R18	0.276 (0.265 - 0.288)	0.869 (0.865 - 0.873)	0.455 (0.442 - 0.468)	0.942 (0.940 - 0.943)	0.352 (0.342 - 0.362)	0.320 (0.309 - 0.331)
	Random Init-R18	0.060 (0.054 - 0.066)	0.726 (0.715 - 0.736)	0.151 (0.137 - 0.166)	0.972 (0.971 - 0.973)	0.119 (0.107 - 0.131)	0.100 (0.088 - 0.112)
	ECGCLIP-R18	0.065 (0.058 - 0.072)	0.736 (0.725 - 0.747)	0.199 (0.181 - 0.216)	0.960 (0.959 - 0.961)	0.125 (0.113 - 0.136)	0.109 (0.096 - 0.121)
	ECGCLIP-R34	0.070 (0.063 - 0.078)	0.759 (0.749 - 0.768)	0.219 (0.202 - 0.236)	0.957 (0.955 - 0.958)	0.129 (0.118 - 0.140)	0.115 (0.104 - 0.127)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AF	Merl-R18	0.057 (0.052 - 0.063)	0.733 (0.722 - 0.743)	0.140 (0.126 - 0.155)	0.972 (0.971 - 0.973)	0.110 (0.099 - 0.122)	0.091 (0.079 - 0.103)
	Random Init-R18	0.830 (0.821 - 0.840)	0.986 (0.985 - 0.987)	0.881 (0.873 - 0.888)	0.977 (0.976 - 0.978)	0.809 (0.803 - 0.816)	0.796 (0.789 - 0.803)
	ECGCLIP-R18	0.860 (0.851 - 0.868)	0.988 (0.987 - 0.989)	0.873 (0.865 - 0.880)	0.980 (0.979 - 0.981)	0.821 (0.814 - 0.827)	0.808 (0.800 - 0.814)
	ECGCLIP-R34	0.900 (0.892 - 0.907)	0.991 (0.990 - 0.992)	0.887 (0.880 - 0.894)	0.985 (0.984 - 0.986)	0.853 (0.847 - 0.858)	0.842 (0.835 - 0.847)
AFL	Merl-R18	0.839 (0.829 - 0.849)	0.987 (0.986 - 0.987)	0.884 (0.877 - 0.891)	0.977 (0.976 - 0.978)	0.811 (0.804 - 0.817)	0.798 (0.791 - 0.804)
	Random Init-R18	0.739 (0.720 - 0.757)	0.987 (0.985 - 0.989)	0.627 (0.607 - 0.649)	0.996 (0.996 - 0.996)	0.675 (0.658 - 0.692)	0.672 (0.655 - 0.690)
	ECGCLIP-R18	0.793 (0.777 - 0.809)	0.991 (0.989 - 0.993)	0.711 (0.691 - 0.730)	0.996 (0.996 - 0.996)	0.733 (0.717 - 0.749)	0.729 (0.712 - 0.745)
	ECGCLIP-R34	0.821 (0.806 - 0.836)	0.992 (0.991 - 0.994)	0.712 (0.692 - 0.733)	0.997 (0.996 - 0.997)	0.748 (0.733 - 0.763)	0.745 (0.730 - 0.760)
JPB	Merl-R18	0.767 (0.749 - 0.784)	0.989 (0.988 - 0.991)	0.649 (0.629 - 0.671)	0.996 (0.996 - 0.997)	0.696 (0.681 - 0.713)	0.693 (0.678 - 0.710)
	Random Init-R18	0.023 (0.020 - 0.027)	0.807 (0.790 - 0.824)	0.165 (0.133 - 0.194)	0.976 (0.975 - 0.977)	0.058 (0.047 - 0.069)	0.066 (0.051 - 0.080)
	ECGCLIP-R18	0.030 (0.025 - 0.035)	0.838 (0.822 - 0.852)	0.128 (0.099 - 0.155)	0.988 (0.987 - 0.988)	0.073 (0.057 - 0.089)	0.073 (0.056 - 0.090)
	ECGCLIP-R34	0.070 (0.059 - 0.083)	0.901 (0.887 - 0.913)	0.254 (0.220 - 0.286)	0.989 (0.989 - 0.990)	0.154 (0.133 - 0.176)	0.161 (0.138 - 0.184)
AJR	Merl-R18	0.023 (0.019 - 0.027)	0.813 (0.797 - 0.829)	0.074 (0.053 - 0.097)	0.990 (0.990 - 0.991)	0.051 (0.036 - 0.066)	0.047 (0.031 - 0.063)
	Random Init-R18	0.145 (0.121 - 0.172)	0.879 (0.867 - 0.891)	0.231 (0.203 - 0.264)	0.993 (0.993 - 0.994)	0.214 (0.189 - 0.242)	0.208 (0.183 - 0.237)
	ECGCLIP-R18	0.153 (0.130 - 0.180)	0.894 (0.883 - 0.905)	0.280 (0.248 - 0.311)	0.992 (0.992 - 0.993)	0.240 (0.213 - 0.267)	0.236 (0.210 - 0.262)
	ECGCLIP-R34	0.172 (0.147 - 0.201)	0.902 (0.892 - 0.913)	0.256 (0.227 - 0.287)	0.995 (0.994 - 0.995)	0.257 (0.229 - 0.287)	0.252 (0.223 - 0.282)
	Merl-R18	0.149 (0.125 - 0.176)	0.882 (0.870 - 0.895)	0.221 (0.193 - 0.252)	0.995 (0.994 - 0.995)	0.225 (0.198 - 0.252)	0.220 (0.192 - 0.247)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
ATach	Random Init-R18	0.093 (0.078 - 0.110)	0.927 (0.919 - 0.935)	0.291 (0.258 - 0.326)	0.986 (0.985 - 0.987)	0.162 (0.142 - 0.180)	0.173 (0.151 - 0.194)
	ECGCLIP-R18	0.153 (0.132 - 0.178)	0.948 (0.941 - 0.954)	0.324 (0.288 - 0.359)	0.991 (0.991 - 0.992)	0.232 (0.206 - 0.257)	0.236 (0.210 - 0.262)
	ECGCLIP-R34	0.222 (0.194 - 0.252)	0.964 (0.957 - 0.969)	0.436 (0.397 - 0.472)	0.992 (0.991 - 0.992)	0.312 (0.283 - 0.336)	0.320 (0.292 - 0.345)
JTach	Merl-R18	0.108 (0.092 - 0.129)	0.934 (0.926 - 0.941)	0.242 (0.210 - 0.273)	0.991 (0.991 - 0.992)	0.178 (0.157 - 0.202)	0.178 (0.156 - 0.202)
	Random Init-R18	0.332 (0.274 - 0.388)	0.973 (0.965 - 0.979)	0.345 (0.293 - 0.399)	0.998 (0.998 - 0.999)	0.365 (0.316 - 0.411)	0.363 (0.315 - 0.410)
	ECGCLIP-R18	0.387 (0.332 - 0.440)	0.980 (0.972 - 0.986)	0.378 (0.325 - 0.434)	0.998 (0.998 - 0.999)	0.396 (0.344 - 0.445)	0.394 (0.344 - 0.444)
JEB	ECGCLIP-R34	0.470 (0.413 - 0.523)	0.984 (0.977 - 0.989)	0.488 (0.432 - 0.542)	0.998 (0.998 - 0.998)	0.472 (0.423 - 0.519)	0.471 (0.422 - 0.518)
	Merl-R18	0.380 (0.324 - 0.435)	0.977 (0.970 - 0.982)	0.415 (0.358 - 0.472)	0.998 (0.997 - 0.998)	0.385 (0.335 - 0.428)	0.384 (0.333 - 0.428)
	Random Init-R18	0.271 (0.232 - 0.320)	0.972 (0.964 - 0.979)	0.493 (0.442 - 0.544)	0.996 (0.996 - 0.997)	0.381 (0.342 - 0.419)	0.389 (0.351 - 0.426)
SVT	ECGCLIP-R18	0.313 (0.268 - 0.370)	0.979 (0.973 - 0.985)	0.493 (0.443 - 0.548)	0.997 (0.997 - 0.998)	0.426 (0.387 - 0.468)	0.428 (0.389 - 0.471)
	ECGCLIP-R34	0.419 (0.367 - 0.478)	0.985 (0.980 - 0.989)	0.422 (0.373 - 0.473)	0.998 (0.998 - 0.999)	0.445 (0.399 - 0.493)	0.444 (0.398 - 0.492)
	Merl-R18	0.267 (0.227 - 0.316)	0.973 (0.966 - 0.980)	0.374 (0.326 - 0.428)	0.998 (0.998 - 0.998)	0.376 (0.332 - 0.422)	0.374 (0.330 - 0.420)
VEB	Random Init-R18	0.631 (0.534 - 0.707)	0.995 (0.989 - 0.999)	0.556 (0.470 - 0.645)	1.000 (0.999 - 1.000)	0.583 (0.507 - 0.651)	0.584 (0.509 - 0.653)
	ECGCLIP-R18	0.677 (0.586 - 0.751)	0.997 (0.993 - 0.999)	0.571 (0.488 - 0.659)	1.000 (1.000 - 1.000)	0.640 (0.560 - 0.710)	0.644 (0.568 - 0.713)
	ECGCLIP-R34	0.757 (0.682 - 0.822)	0.994 (0.982 - 1.000)	0.627 (0.540 - 0.706)	1.000 (1.000 - 1.000)	0.702 (0.629 - 0.766)	0.707 (0.638 - 0.767)
	Merl-R18	0.651 (0.563 - 0.724)	0.996 (0.991 - 0.999)	0.524 (0.439 - 0.610)	1.000 (1.000 - 1.000)	0.595 (0.513 - 0.670)	0.600 (0.522 - 0.674)
	Random Init-R18	0.100 (0.069 - 0.141)	0.966 (0.955 - 0.976)	0.205 (0.144 - 0.267)	0.999 (0.999 - 0.999)	0.202 (0.140 - 0.257)	0.201 (0.139 - 0.256)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
VTach	ECGCLIP-R18	0.146 (0.101 - 0.205)	0.971 (0.962 - 0.980)	0.224 (0.160 - 0.288)	0.999 (0.999 - 0.999)	0.260 (0.194 - 0.327)	0.262 (0.195 - 0.330)
	ECGCLIP-R34	0.201 (0.138 - 0.266)	0.974 (0.965 - 0.982)	0.217 (0.156 - 0.279)	0.999 (0.999 - 1.000)	0.270 (0.198 - 0.336)	0.278 (0.206 - 0.342)
	Merl-R18	0.092 (0.065 - 0.126)	0.968 (0.956 - 0.977)	0.391 (0.320 - 0.469)	0.995 (0.995 - 0.995)	0.167 (0.132 - 0.202)	0.202 (0.162 - 0.242)
	Random Init-R18	0.270 (0.193 - 0.356)	0.975 (0.966 - 0.982)	0.396 (0.316 - 0.477)	0.999 (0.999 - 0.999)	0.372 (0.302 - 0.441)	0.372 (0.302 - 0.441)
	ECGCLIP-R18	0.296 (0.215 - 0.384)	0.979 (0.972 - 0.984)	0.389 (0.310 - 0.467)	0.999 (0.999 - 0.999)	0.372 (0.299 - 0.442)	0.371 (0.299 - 0.443)
	ECGCLIP-R34	0.308 (0.228 - 0.400)	0.982 (0.976 - 0.986)	0.389 (0.309 - 0.467)	0.999 (0.999 - 0.999)	0.368 (0.299 - 0.435)	0.368 (0.299 - 0.435)
	Merl-R18	0.276 (0.201 - 0.363)	0.975 (0.967 - 0.982)	0.242 (0.177 - 0.311)	1.000 (1.000 - 1.000)	0.333 (0.249 - 0.411)	0.360 (0.275 - 0.438)
	Random Init-R18	0.955 (0.951 - 0.958)	0.996 (0.995 - 0.996)	0.878 (0.869 - 0.885)	0.992 (0.992 - 0.993)	0.881 (0.875 - 0.887)	0.873 (0.867 - 0.879)
	ECGCLIP-R18	0.960 (0.957 - 0.963)	0.996 (0.996 - 0.997)	0.879 (0.871 - 0.886)	0.993 (0.993 - 0.994)	0.888 (0.883 - 0.894)	0.881 (0.875 - 0.886)
	ECGCLIP-R34	0.963 (0.959 - 0.965)	0.997 (0.996 - 0.997)	0.887 (0.879 - 0.894)	0.993 (0.993 - 0.994)	0.892 (0.887 - 0.897)	0.885 (0.879 - 0.890)
RBBB	Merl-R18	0.957 (0.953 - 0.960)	0.996 (0.996 - 0.997)	0.887 (0.879 - 0.894)	0.992 (0.991 - 0.992)	0.884 (0.878 - 0.889)	0.876 (0.870 - 0.882)
	Random Init-R18	0.850 (0.840 - 0.860)	0.990 (0.989 - 0.991)	0.787 (0.774 - 0.800)	0.991 (0.991 - 0.992)	0.778 (0.768 - 0.789)	0.770 (0.760 - 0.781)
	ECGCLIP-R18	0.863 (0.853 - 0.872)	0.992 (0.990 - 0.993)	0.810 (0.798 - 0.822)	0.991 (0.990 - 0.992)	0.788 (0.778 - 0.797)	0.780 (0.770 - 0.790)
	ECGCLIP-R34	0.878 (0.868 - 0.887)	0.992 (0.991 - 0.993)	0.792 (0.779 - 0.805)	0.993 (0.993 - 0.994)	0.800 (0.789 - 0.810)	0.793 (0.782 - 0.803)
1° AVB	Merl-R18	0.855 (0.845 - 0.864)	0.990 (0.989 - 0.991)	0.823 (0.810 - 0.835)	0.989 (0.988 - 0.990)	0.776 (0.766 - 0.786)	0.769 (0.759 - 0.779)
	Random Init-R18	0.227 (0.204 - 0.250)	0.932 (0.926 - 0.938)	0.298 (0.273 - 0.324)	0.992 (0.991 - 0.992)	0.293 (0.270 - 0.315)	0.285 (0.262 - 0.307)
	ECGCLIP-R18	0.261 (0.236 - 0.286)	0.939 (0.934 - 0.945)	0.353 (0.327 - 0.379)	0.991 (0.991 - 0.992)	0.330 (0.309 - 0.351)	0.323 (0.301 - 0.344)
QT Prolong							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
ER	ECGCLIP-R34	0.299 (0.273 - 0.324)	0.947 (0.943 - 0.952)	0.448 (0.422 - 0.475)	0.987 (0.987 - 0.988)	0.346 (0.326 - 0.366)	0.347 (0.326 - 0.366)
	Merl-R18	0.240 (0.216 - 0.263)	0.933 (0.927 - 0.938)	0.320 (0.294 - 0.345)	0.991 (0.990 - 0.991)	0.301 (0.279 - 0.321)	0.293 (0.271 - 0.313)
	Random Init-R18	0.230 (0.188 - 0.270)	0.952 (0.943 - 0.960)	0.309 (0.268 - 0.350)	0.997 (0.997 - 0.997)	0.308 (0.269 - 0.346)	0.305 (0.266 - 0.344)
	ECGCLIP-R18	0.258 (0.215 - 0.303)	0.964 (0.957 - 0.971)	0.255 (0.215 - 0.295)	0.998 (0.998 - 0.999)	0.310 (0.266 - 0.352)	0.316 (0.271 - 0.356)
	ECGCLIP-R34	0.262 (0.222 - 0.306)	0.968 (0.961 - 0.974)	0.383 (0.338 - 0.430)	0.997 (0.996 - 0.997)	0.351 (0.312 - 0.388)	0.349 (0.310 - 0.387)
	Merl-R18	0.225 (0.184 - 0.266)	0.941 (0.929 - 0.951)	0.338 (0.294 - 0.384)	0.996 (0.995 - 0.996)	0.283 (0.246 - 0.320)	0.283 (0.246 - 0.322)
	Random Init-R18	0.443 (0.405 - 0.480)	0.986 (0.984 - 0.988)	0.577 (0.545 - 0.612)	0.994 (0.994 - 0.995)	0.492 (0.464 - 0.520)	0.493 (0.465 - 0.521)
	ECGCLIP-R18	0.513 (0.474 - 0.551)	0.991 (0.989 - 0.992)	0.577 (0.543 - 0.611)	0.996 (0.995 - 0.996)	0.542 (0.514 - 0.571)	0.540 (0.511 - 0.569)
	ECGCLIP-R34	0.538 (0.499 - 0.576)	0.992 (0.991 - 0.993)	0.623 (0.588 - 0.655)	0.995 (0.995 - 0.996)	0.554 (0.526 - 0.580)	0.553 (0.527 - 0.580)
LAFB	Merl-R18	0.482 (0.443 - 0.518)	0.989 (0.987 - 0.991)	0.545 (0.513 - 0.581)	0.995 (0.995 - 0.996)	0.505 (0.477 - 0.534)	0.503 (0.474 - 0.531)
	Random Init-R18	0.921 (0.905 - 0.935)	0.999 (0.999 - 0.999)	0.851 (0.829 - 0.875)	0.999 (0.999 - 0.999)	0.856 (0.837 - 0.873)	0.854 (0.836 - 0.872)
	ECGCLIP-R18	0.929 (0.914 - 0.941)	0.999 (0.998 - 0.999)	0.842 (0.820 - 0.865)	0.999 (0.999 - 0.999)	0.855 (0.837 - 0.873)	0.854 (0.836 - 0.872)
	ECGCLIP-R34	0.935 (0.922 - 0.947)	0.999 (0.998 - 0.999)	0.832 (0.808 - 0.857)	0.999 (0.999 - 0.999)	0.854 (0.835 - 0.872)	0.853 (0.834 - 0.871)
LBBB	Merl-R18	0.913 (0.896 - 0.928)	0.999 (0.999 - 0.999)	0.825 (0.799 - 0.849)	0.999 (0.999 - 0.999)	0.842 (0.822 - 0.860)	0.841 (0.821 - 0.859)
	Random Init-R18	0.624 (0.588 - 0.661)	0.992 (0.991 - 0.994)	0.635 (0.600 - 0.668)	0.997 (0.996 - 0.997)	0.614 (0.586 - 0.640)	0.611 (0.583 - 0.638)
	ECGCLIP-R18	0.661 (0.626 - 0.696)	0.994 (0.992 - 0.995)	0.635 (0.600 - 0.668)	0.997 (0.997 - 0.997)	0.635 (0.607 - 0.662)	0.632 (0.604 - 0.659)
IVB	ECGCLIP-R34	0.691 (0.657 - 0.726)	0.995 (0.994 - 0.996)	0.629 (0.594 - 0.662)	0.998 (0.997 - 0.998)	0.654 (0.626 - 0.680)	0.652 (0.624 - 0.679)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
Short PR	Merl-R18	0.618 (0.581 - 0.656)	0.992 (0.991 - 0.994)	0.620 (0.587 - 0.655)	0.997 (0.997 - 0.997)	0.623 (0.594 - 0.653)	0.620 (0.591 - 0.650)
	Random Init-R18	0.211 (0.164 - 0.268)	0.986 (0.977 - 0.992)	0.375 (0.306 - 0.446)	0.998 (0.998 - 0.998)	0.300 (0.246 - 0.352)	0.305 (0.249 - 0.358)
	ECGCLIP-R18	0.244 (0.190 - 0.305)	0.988 (0.980 - 0.993)	0.370 (0.302 - 0.440)	0.998 (0.998 - 0.998)	0.289 (0.236 - 0.344)	0.295 (0.242 - 0.350)
	ECGCLIP-R34	0.269 (0.216 - 0.335)	0.990 (0.985 - 0.994)	0.458 (0.389 - 0.528)	0.998 (0.998 - 0.998)	0.368 (0.312 - 0.424)	0.374 (0.318 - 0.429)
VPE	Merl-R18	0.206 (0.157 - 0.261)	0.982 (0.969 - 0.991)	0.490 (0.420 - 0.559)	0.997 (0.997 - 0.997)	0.301 (0.255 - 0.349)	0.325 (0.278 - 0.371)
	Random Init-R18	0.789 (0.732 - 0.835)	0.985 (0.973 - 0.995)	0.639 (0.571 - 0.697)	1.000 (1.000 - 1.000)	0.739 (0.682 - 0.786)	0.748 (0.694 - 0.793)
	ECGCLIP-R18	0.826 (0.775 - 0.871)	0.989 (0.980 - 0.996)	0.703 (0.636 - 0.761)	1.000 (1.000 - 1.000)	0.787 (0.733 - 0.830)	0.792 (0.741 - 0.834)
	ECGCLIP-R34	0.823 (0.772 - 0.869)	0.993 (0.988 - 0.997)	0.683 (0.615 - 0.742)	1.000 (1.000 - 1.000)	0.795 (0.745 - 0.837)	0.806 (0.761 - 0.843)
3° AVB	Merl-R18	0.759 (0.697 - 0.808)	0.987 (0.979 - 0.994)	0.693 (0.626 - 0.755)	1.000 (1.000 - 1.000)	0.747 (0.693 - 0.794)	0.748 (0.695 - 0.795)
	Random Init-R18	0.157 (0.127 - 0.195)	0.987 (0.984 - 0.990)	0.365 (0.309 - 0.427)	0.996 (0.996 - 0.996)	0.240 (0.199 - 0.279)	0.253 (0.211 - 0.294)
	ECGCLIP-R18	0.228 (0.181 - 0.281)	0.990 (0.987 - 0.992)	0.550 (0.486 - 0.611)	0.995 (0.994 - 0.995)	0.293 (0.253 - 0.331)	0.329 (0.290 - 0.369)
	ECGCLIP-R34	0.347 (0.282 - 0.411)	0.991 (0.987 - 0.994)	0.289 (0.238 - 0.344)	0.999 (0.999 - 0.999)	0.336 (0.280 - 0.392)	0.340 (0.283 - 0.395)
2° 1 Type AVB	Merl-R18	0.145 (0.116 - 0.178)	0.988 (0.985 - 0.990)	0.373 (0.311 - 0.431)	0.996 (0.996 - 0.996)	0.239 (0.197 - 0.278)	0.254 (0.211 - 0.294)
	Random Init-R18	0.125 (0.092 - 0.168)	0.987 (0.981 - 0.992)	0.154 (0.099 - 0.215)	0.999 (0.999 - 0.999)	0.170 (0.109 - 0.230)	0.170 (0.109 - 0.232)
	ECGCLIP-R18	0.156 (0.119 - 0.202)	0.993 (0.991 - 0.995)	0.412 (0.329 - 0.493)	0.998 (0.998 - 0.998)	0.283 (0.225 - 0.340)	0.297 (0.238 - 0.354)
	ECGCLIP-R34	0.204 (0.162 - 0.259)	0.995 (0.994 - 0.996)	0.434 (0.351 - 0.519)	0.999 (0.998 - 0.999)	0.342 (0.281 - 0.404)	0.349 (0.286 - 0.412)
	Merl-R18	0.116 (0.084 - 0.158)	0.988 (0.982 - 0.992)	0.147 (0.089 - 0.208)	0.999 (0.999 - 0.999)	0.170 (0.106 - 0.237)	0.171 (0.106 - 0.239)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
2° 2 Type AVB	Random Init-R18	0.084 (0.051 - 0.137)	0.994 (0.990 - 0.997)	0.240 (0.128 - 0.365)	0.999 (0.999 - 0.999)	0.170 (0.087 - 0.256)	0.177 (0.094 - 0.267)
	ECGCLIP-R18	0.232 (0.145 - 0.354)	0.997 (0.994 - 0.999)	0.420 (0.277 - 0.556)	1.000 (0.999 - 1.000)	0.362 (0.242 - 0.478)	0.365 (0.244 - 0.479)
	ECGCLIP-R34	0.315 (0.220 - 0.435)	0.999 (0.998 - 1.000)	0.500 (0.364 - 0.639)	1.000 (0.999 - 1.000)	0.410 (0.296 - 0.515)	0.416 (0.303 - 0.524)
	Merl-R18	0.118 (0.067 - 0.188)	0.994 (0.991 - 0.997)	0.220 (0.113 - 0.339)	1.000 (1.000 - 1.000)	0.232 (0.121 - 0.337)	0.232 (0.122 - 0.337)
LVH	Random Init-R18	0.884 (0.877 - 0.891)	0.990 (0.989 - 0.991)	0.798 (0.788 - 0.808)	0.989 (0.988 - 0.989)	0.805 (0.797 - 0.813)	0.793 (0.785 - 0.801)
	ECGCLIP-R18	0.895 (0.888 - 0.902)	0.992 (0.991 - 0.992)	0.800 (0.790 - 0.810)	0.990 (0.990 - 0.991)	0.816 (0.809 - 0.823)	0.805 (0.798 - 0.813)
	ECGCLIP-R34	0.890 (0.883 - 0.896)	0.991 (0.990 - 0.992)	0.810 (0.800 - 0.820)	0.989 (0.988 - 0.990)	0.813 (0.806 - 0.821)	0.802 (0.794 - 0.810)
	Merl-R18	0.878 (0.871 - 0.885)	0.989 (0.988 - 0.990)	0.793 (0.784 - 0.804)	0.988 (0.988 - 0.989)	0.800 (0.792 - 0.808)	0.788 (0.780 - 0.796)
LAH	Random Init-R18	0.498 (0.475 - 0.524)	0.974 (0.971 - 0.976)	0.574 (0.553 - 0.595)	0.988 (0.987 - 0.989)	0.510 (0.492 - 0.529)	0.504 (0.486 - 0.522)
	ECGCLIP-R18	0.516 (0.493 - 0.541)	0.977 (0.975 - 0.979)	0.576 (0.555 - 0.601)	0.989 (0.989 - 0.990)	0.527 (0.509 - 0.548)	0.520 (0.502 - 0.541)
	ECGCLIP-R34	0.528 (0.506 - 0.554)	0.978 (0.976 - 0.980)	0.644 (0.622 - 0.665)	0.985 (0.985 - 0.986)	0.520 (0.503 - 0.538)	0.520 (0.503 - 0.538)
	Merl-R18	0.497 (0.474 - 0.524)	0.975 (0.973 - 0.977)	0.607 (0.586 - 0.630)	0.986 (0.986 - 0.987)	0.511 (0.494 - 0.530)	0.508 (0.491 - 0.527)
RVH	Random Init-R18	0.439 (0.407 - 0.472)	0.981 (0.978 - 0.984)	0.513 (0.484 - 0.544)	0.993 (0.993 - 0.994)	0.453 (0.428 - 0.480)	0.450 (0.426 - 0.478)
	ECGCLIP-R18	0.496 (0.464 - 0.529)	0.986 (0.983 - 0.989)	0.555 (0.524 - 0.585)	0.994 (0.994 - 0.995)	0.511 (0.485 - 0.537)	0.508 (0.481 - 0.534)
	ECGCLIP-R34	0.538 (0.504 - 0.572)	0.988 (0.985 - 0.990)	0.593 (0.564 - 0.623)	0.995 (0.994 - 0.995)	0.542 (0.516 - 0.568)	0.540 (0.513 - 0.566)
	Merl-R18	0.450 (0.418 - 0.485)	0.983 (0.980 - 0.985)	0.457 (0.428 - 0.487)	0.995 (0.995 - 0.996)	0.463 (0.435 - 0.490)	0.458 (0.430 - 0.486)
VA	Random Init-R18	0.202 (0.158 - 0.258)	0.987 (0.982 - 0.991)	0.471 (0.399 - 0.540)	0.997 (0.997 - 0.998)	0.302 (0.250 - 0.348)	0.322 (0.270 - 0.368)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RAH	ECGCLIP-R18	0.227 (0.171 - 0.287)	0.990 (0.987 - 0.993)	0.431 (0.358 - 0.503)	0.998 (0.997 - 0.998)	0.301 (0.248 - 0.350)	0.314 (0.259 - 0.366)
	ECGCLIP-R34	0.252 (0.195 - 0.314)	0.992 (0.989 - 0.994)	0.420 (0.352 - 0.491)	0.998 (0.998 - 0.998)	0.326 (0.273 - 0.381)	0.333 (0.281 - 0.388)
	Merl-R18	0.198 (0.150 - 0.256)	0.989 (0.985 - 0.992)	0.443 (0.372 - 0.511)	0.997 (0.997 - 0.997)	0.273 (0.227 - 0.321)	0.294 (0.247 - 0.342)
	Random Init-R18	0.024 (0.015 - 0.045)	0.941 (0.920 - 0.960)	0.090 (0.035 - 0.151)	0.999 (0.998 - 0.999)	0.067 (0.026 - 0.113)	0.068 (0.026 - 0.115)
	ECGCLIP-R18	0.024 (0.015 - 0.041)	0.957 (0.943 - 0.969)	0.101 (0.043 - 0.172)	0.999 (0.999 - 0.999)	0.079 (0.035 - 0.130)	0.080 (0.034 - 0.133)
	ECGCLIP-R34	0.038 (0.023 - 0.063)	0.967 (0.952 - 0.978)	0.056 (0.012 - 0.108)	1.000 (0.999 - 1.000)	0.069 (0.014 - 0.130)	0.070 (0.014 - 0.136)
	Merl-R18	0.022 (0.013 - 0.048)	0.938 (0.917 - 0.957)	0.011 (0.000 - 0.037)	1.000 (1.000 - 1.000)	0.015 (0.000 - 0.048)	0.016 (-0.001 - 0.051)
	Random Init-R18	0.000 (0.000 - 0.001)	0.715 (0.613 - 0.807)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.001 (0.000 - 0.001)	0.705 (0.581 - 0.808)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.001 (0.000 - 0.001)	0.741 (0.650 - 0.825)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
DC	Merl-R18	0.000 (0.000 - 0.001)	0.694 (0.588 - 0.788)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	Random Init-R18	0.022 (0.012 - 0.041)	0.977 (0.966 - 0.986)	0.068 (0.000 - 0.161)	0.999 (0.999 - 0.999)	0.042 (0.000 - 0.096)	0.045 (-0.001 - 0.103)
	ECGCLIP-R18	0.022 (0.011 - 0.047)	0.978 (0.970 - 0.986)	0.091 (0.018 - 0.195)	0.999 (0.999 - 0.999)	0.051 (0.012 - 0.107)	0.056 (0.012 - 0.119)
	ECGCLIP-R34	0.044 (0.019 - 0.098)	0.982 (0.974 - 0.989)	0.114 (0.024 - 0.213)	1.000 (1.000 - 1.000)	0.143 (0.031 - 0.258)	0.148 (0.032 - 0.272)
BVH	Merl-R18	0.016 (0.008 - 0.036)	0.967 (0.949 - 0.981)	0.045 (0.000 - 0.122)	1.000 (1.000 - 1.000)	0.056 (0.000 - 0.146)	0.058 (0.000 - 0.151)
	Random Init-R18	0.643 (0.606 - 0.678)	0.993 (0.992 - 0.994)	0.751 (0.722 - 0.778)	0.996 (0.995 - 0.996)	0.660 (0.636 - 0.682)	0.662 (0.638 - 0.684)
	ECGCLIP-R18	0.686 (0.651 - 0.719)	0.995 (0.993 - 0.996)	0.773 (0.745 - 0.797)	0.996 (0.996 - 0.997)	0.694 (0.671 - 0.715)	0.695 (0.672 - 0.716)
VVI							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
VAT	ECGCLIP-R34	0.766 (0.732 - 0.796)	0.996 (0.995 - 0.997)	0.834 (0.809 - 0.857)	0.997 (0.997 - 0.997)	0.764 (0.742 - 0.784)	0.764 (0.744 - 0.784)
	Merl-R18	0.675 (0.640 - 0.710)	0.994 (0.993 - 0.995)	0.738 (0.711 - 0.768)	0.997 (0.997 - 0.997)	0.695 (0.671 - 0.719)	0.694 (0.669 - 0.718)
	Random Init-R18	0.644 (0.597 - 0.692)	0.990 (0.985 - 0.994)	0.582 (0.536 - 0.631)	0.999 (0.999 - 0.999)	0.622 (0.580 - 0.662)	0.622 (0.581 - 0.663)
	ECGCLIP-R18	0.742 (0.698 - 0.781)	0.994 (0.992 - 0.996)	0.657 (0.609 - 0.700)	0.999 (0.999 - 0.999)	0.706 (0.667 - 0.742)	0.707 (0.668 - 0.743)
	ECGCLIP-R34	0.806 (0.770 - 0.841)	0.994 (0.991 - 0.996)	0.698 (0.653 - 0.739)	1.000 (0.999 - 1.000)	0.771 (0.737 - 0.803)	0.774 (0.740 - 0.806)
	Merl-R18	0.707 (0.660 - 0.753)	0.995 (0.993 - 0.996)	0.675 (0.628 - 0.719)	0.999 (0.999 - 0.999)	0.683 (0.643 - 0.718)	0.681 (0.642 - 0.717)
DDD	Random Init-R18	0.357 (0.292 - 0.434)	0.992 (0.989 - 0.994)	0.402 (0.340 - 0.466)	0.999 (0.998 - 0.999)	0.392 (0.335 - 0.448)	0.391 (0.334 - 0.447)
	ECGCLIP-R18	0.492 (0.427 - 0.562)	0.994 (0.993 - 0.996)	0.511 (0.451 - 0.579)	0.999 (0.998 - 0.999)	0.476 (0.423 - 0.530)	0.476 (0.423 - 0.531)
	ECGCLIP-R34	0.666 (0.608 - 0.727)	0.996 (0.995 - 0.998)	0.589 (0.527 - 0.655)	0.999 (0.999 - 1.000)	0.620 (0.566 - 0.673)	0.620 (0.565 - 0.674)
	Merl-R18	0.399 (0.330 - 0.470)	0.994 (0.992 - 0.995)	0.438 (0.377 - 0.500)	0.998 (0.998 - 0.999)	0.388 (0.336 - 0.444)	0.389 (0.337 - 0.445)
AAI	Random Init-R18	0.190 (0.146 - 0.251)	0.975 (0.964 - 0.983)	0.255 (0.200 - 0.321)	0.999 (0.999 - 0.999)	0.261 (0.206 - 0.320)	0.260 (0.205 - 0.319)
	ECGCLIP-R18	0.407 (0.334 - 0.491)	0.982 (0.970 - 0.991)	0.281 (0.218 - 0.345)	1.000 (1.000 - 1.000)	0.379 (0.303 - 0.447)	0.404 (0.329 - 0.471)
	ECGCLIP-R34	0.705 (0.639 - 0.764)	0.990 (0.978 - 0.997)	0.541 (0.467 - 0.609)	1.000 (1.000 - 1.000)	0.639 (0.572 - 0.697)	0.649 (0.585 - 0.704)
	Merl-R18	0.360 (0.290 - 0.439)	0.984 (0.972 - 0.992)	0.459 (0.394 - 0.530)	0.999 (0.998 - 0.999)	0.407 (0.351 - 0.465)	0.409 (0.352 - 0.466)
VP	Random Init-R18	0.098 (0.059 - 0.159)	0.979 (0.972 - 0.985)	0.538 (0.446 - 0.628)	0.992 (0.992 - 0.993)	0.123 (0.095 - 0.150)	0.191 (0.155 - 0.226)
	ECGCLIP-R18	0.138 (0.085 - 0.211)	0.981 (0.975 - 0.987)	0.513 (0.423 - 0.600)	0.994 (0.993 - 0.994)	0.140 (0.109 - 0.170)	0.202 (0.163 - 0.240)
	ECGCLIP-R34	0.117 (0.072 - 0.179)	0.983 (0.976 - 0.989)	0.291 (0.212 - 0.372)	0.998 (0.997 - 0.998)	0.165 (0.116 - 0.215)	0.182 (0.129 - 0.235)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
PSM	Merl-R18	0.102 (0.065 - 0.160)	0.983 (0.978 - 0.988)	0.538 (0.449 - 0.630)	0.992 (0.991 - 0.992)	0.116 (0.091 - 0.143)	0.185 (0.151 - 0.220)
	Random Init-R18	0.024 (0.011 - 0.056)	0.985 (0.978 - 0.990)	0.057 (0.000 - 0.152)	1.000 (1.000 - 1.000)	0.075 (0.000 - 0.191)	0.079 (0.000 - 0.203)
	ECGCLIP-R18	0.031 (0.012 - 0.102)	0.984 (0.976 - 0.991)	0.057 (0.000 - 0.152)	1.000 (1.000 - 1.000)	0.078 (0.000 - 0.196)	0.084 (0.000 - 0.207)
	ECGCLIP-R34	0.069 (0.033 - 0.147)	0.991 (0.986 - 0.996)	0.086 (0.000 - 0.200)	1.000 (1.000 - 1.000)	0.100 (0.000 - 0.219)	0.101 (0.000 - 0.220)
HK	Merl-R18	0.024 (0.012 - 0.055)	0.988 (0.981 - 0.992)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Random Init-R18	0.003 (0.002 - 0.007)	0.865 (0.810 - 0.912)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.001 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.005 (0.002 - 0.011)	0.852 (0.794 - 0.908)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.001 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.006 (0.003 - 0.013)	0.907 (0.872 - 0.940)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.001 - 0.001)	0.000 (0.000 - 0.000)
LK	Merl-R18	0.003 (0.002 - 0.005)	0.849 (0.796 - 0.893)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.001 - 0.001)	0.000 (0.000 - 0.000)
	Random Init-R18	0.005 (0.003 - 0.010)	0.911 (0.880 - 0.939)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.001 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.007 (0.004 - 0.014)	0.920 (0.888 - 0.949)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.001 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.011 (0.006 - 0.022)	0.926 (0.887 - 0.958)	0.018 (0.000 - 0.064)	1.000 (0.999 - 1.000)	0.019 (0.000 - 0.067)	0.019 (-0.001 - 0.066)
	Merl-R18	0.005 (0.003 - 0.008)	0.891 (0.851 - 0.930)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.001 - 0.001)	0.000 (0.000 - 0.000)
Tier 2: ECHO							
Diagnosis							
NER	Random Init-R18	0.858 (0.853 - 0.863)	0.893 (0.890 - 0.896)	0.885 (0.881 - 0.889)	0.725 (0.720 - 0.730)	0.797 (0.793 - 0.801)	0.609 (0.603 - 0.616)
	ECGCLIP-R18	0.871 (0.866 - 0.875)	0.902 (0.900 - 0.905)	0.887 (0.883 - 0.892)	0.743 (0.737 - 0.748)	0.806 (0.802 - 0.810)	0.629 (0.622 - 0.635)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAE	ECGCLIP-R34	0.875 (0.870 - 0.879)	0.905 (0.902 - 0.907)	0.894 (0.890 - 0.898)	0.739 (0.734 - 0.744)	0.808 (0.804 - 0.812)	0.633 (0.626 - 0.639)
	Merl-R18	0.866 (0.861 - 0.870)	0.899 (0.896 - 0.901)	0.885 (0.881 - 0.889)	0.740 (0.735 - 0.745)	0.803 (0.800 - 0.807)	0.624 (0.617 - 0.630)
	Random Init-R18	0.425 (0.414 - 0.436)	0.798 (0.794 - 0.803)	0.651 (0.641 - 0.662)	0.770 (0.765 - 0.774)	0.469 (0.461 - 0.477)	0.344 (0.335 - 0.354)
	ECGCLIP-R18	0.455 (0.444 - 0.466)	0.813 (0.809 - 0.818)	0.614 (0.604 - 0.626)	0.816 (0.813 - 0.820)	0.490 (0.481 - 0.498)	0.371 (0.361 - 0.381)
	ECGCLIP-R34	0.464 (0.452 - 0.476)	0.817 (0.813 - 0.822)	0.627 (0.617 - 0.637)	0.815 (0.811 - 0.819)	0.496 (0.487 - 0.504)	0.378 (0.369 - 0.388)
	Merl-R18	0.442 (0.430 - 0.453)	0.807 (0.802 - 0.812)	0.611 (0.601 - 0.622)	0.810 (0.806 - 0.814)	0.482 (0.473 - 0.490)	0.360 (0.351 - 0.370)
HCM	Random Init-R18	0.455 (0.442 - 0.469)	0.834 (0.829 - 0.840)	0.545 (0.533 - 0.558)	0.881 (0.878 - 0.884)	0.448 (0.438 - 0.458)	0.367 (0.355 - 0.378)
	ECGCLIP-R18	0.475 (0.462 - 0.489)	0.845 (0.840 - 0.850)	0.533 (0.520 - 0.546)	0.897 (0.895 - 0.900)	0.463 (0.453 - 0.474)	0.385 (0.374 - 0.397)
	ECGCLIP-R34	0.484 (0.470 - 0.497)	0.847 (0.842 - 0.852)	0.552 (0.540 - 0.565)	0.894 (0.891 - 0.897)	0.471 (0.460 - 0.481)	0.393 (0.382 - 0.405)
	Merl-R18	0.462 (0.449 - 0.476)	0.838 (0.833 - 0.843)	0.502 (0.490 - 0.515)	0.905 (0.902 - 0.907)	0.454 (0.443 - 0.464)	0.375 (0.363 - 0.386)
nHCM	Random Init-R18	0.363 (0.351 - 0.377)	0.820 (0.815 - 0.826)	0.528 (0.515 - 0.542)	0.869 (0.866 - 0.872)	0.407 (0.397 - 0.417)	0.326 (0.315 - 0.338)
	ECGCLIP-R18	0.379 (0.366 - 0.394)	0.830 (0.825 - 0.836)	0.549 (0.536 - 0.564)	0.871 (0.868 - 0.875)	0.423 (0.412 - 0.434)	0.346 (0.334 - 0.357)
	ECGCLIP-R34	0.386 (0.373 - 0.400)	0.833 (0.827 - 0.838)	0.589 (0.576 - 0.604)	0.858 (0.854 - 0.861)	0.429 (0.419 - 0.439)	0.355 (0.343 - 0.366)
BAE	Merl-R18	0.368 (0.356 - 0.383)	0.824 (0.818 - 0.829)	0.523 (0.509 - 0.537)	0.880 (0.877 - 0.883)	0.418 (0.408 - 0.429)	0.340 (0.328 - 0.351)
	Random Init-R18	0.414 (0.392 - 0.437)	0.926 (0.921 - 0.931)	0.549 (0.527 - 0.569)	0.963 (0.961 - 0.965)	0.472 (0.455 - 0.489)	0.449 (0.430 - 0.465)
	ECGCLIP-R18	0.453 (0.431 - 0.477)	0.933 (0.929 - 0.938)	0.573 (0.551 - 0.593)	0.965 (0.963 - 0.967)	0.496 (0.479 - 0.513)	0.473 (0.456 - 0.491)
	ECGCLIP-R34	0.466 (0.443 - 0.491)	0.935 (0.931 - 0.940)	0.546 (0.524 - 0.568)	0.969 (0.968 - 0.971)	0.498 (0.480 - 0.515)	0.474 (0.456 - 0.492)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
HFrEF	Merl-R18	0.437 (0.415 - 0.461)	0.930 (0.925 - 0.935)	0.576 (0.554 - 0.597)	0.961 (0.960 - 0.963)	0.482 (0.465 - 0.499)	0.460 (0.443 - 0.478)
	Random Init-R18	0.506 (0.477 - 0.538)	0.956 (0.950 - 0.961)	0.744 (0.721 - 0.768)	0.964 (0.962 - 0.966)	0.525 (0.505 - 0.544)	0.530 (0.511 - 0.549)
	ECGCLIP-R18	0.586 (0.553 - 0.615)	0.962 (0.957 - 0.968)	0.705 (0.680 - 0.727)	0.975 (0.974 - 0.977)	0.578 (0.557 - 0.598)	0.571 (0.550 - 0.591)
	ECGCLIP-R34	0.584 (0.552 - 0.616)	0.964 (0.958 - 0.968)	0.745 (0.720 - 0.769)	0.972 (0.971 - 0.974)	0.579 (0.558 - 0.598)	0.577 (0.557 - 0.596)
RMWSSFLV	Merl-R18	0.556 (0.524 - 0.587)	0.959 (0.953 - 0.964)	0.709 (0.683 - 0.734)	0.972 (0.970 - 0.974)	0.557 (0.538 - 0.578)	0.553 (0.533 - 0.574)
	Random Init-R18	0.453 (0.425 - 0.484)	0.947 (0.941 - 0.953)	0.466 (0.439 - 0.493)	0.988 (0.988 - 0.990)	0.494 (0.470 - 0.517)	0.482 (0.458 - 0.506)
	ECGCLIP-R18	0.489 (0.458 - 0.519)	0.956 (0.950 - 0.960)	0.539 (0.512 - 0.566)	0.986 (0.985 - 0.987)	0.521 (0.498 - 0.543)	0.508 (0.484 - 0.531)
	ECGCLIP-R34	0.502 (0.472 - 0.531)	0.958 (0.953 - 0.963)	0.549 (0.521 - 0.576)	0.985 (0.984 - 0.987)	0.527 (0.504 - 0.549)	0.514 (0.491 - 0.537)
DCM	Merl-R18	0.479 (0.450 - 0.509)	0.953 (0.948 - 0.958)	0.496 (0.468 - 0.522)	0.988 (0.987 - 0.989)	0.509 (0.485 - 0.532)	0.497 (0.472 - 0.520)
	Random Init-R18	0.376 (0.337 - 0.415)	0.968 (0.961 - 0.974)	0.468 (0.429 - 0.510)	0.990 (0.990 - 0.991)	0.418 (0.384 - 0.450)	0.412 (0.378 - 0.445)
	ECGCLIP-R18	0.389 (0.348 - 0.431)	0.971 (0.964 - 0.977)	0.622 (0.585 - 0.659)	0.985 (0.984 - 0.986)	0.445 (0.415 - 0.475)	0.455 (0.427 - 0.485)
	ECGCLIP-R34	0.392 (0.353 - 0.434)	0.974 (0.967 - 0.979)	0.560 (0.518 - 0.601)	0.988 (0.987 - 0.989)	0.447 (0.416 - 0.476)	0.448 (0.418 - 0.477)
oHCM	Merl-R18	0.368 (0.329 - 0.406)	0.969 (0.961 - 0.975)	0.470 (0.430 - 0.510)	0.991 (0.990 - 0.992)	0.426 (0.390 - 0.458)	0.420 (0.384 - 0.453)
	Random Init-R18	0.250 (0.210 - 0.293)	0.947 (0.935 - 0.957)	0.358 (0.309 - 0.405)	0.992 (0.991 - 0.993)	0.312 (0.273 - 0.350)	0.308 (0.268 - 0.348)
	ECGCLIP-R18	0.292 (0.249 - 0.339)	0.957 (0.946 - 0.966)	0.363 (0.315 - 0.412)	0.994 (0.994 - 0.995)	0.364 (0.320 - 0.405)	0.359 (0.315 - 0.400)
	ECGCLIP-R34	0.319 (0.276 - 0.368)	0.959 (0.947 - 0.968)	0.495 (0.444 - 0.544)	0.992 (0.991 - 0.993)	0.407 (0.364 - 0.448)	0.407 (0.364 - 0.449)
	Merl-R18	0.264 (0.224 - 0.309)	0.950 (0.939 - 0.960)	0.414 (0.368 - 0.461)	0.992 (0.991 - 0.993)	0.355 (0.316 - 0.392)	0.352 (0.313 - 0.390)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
ICM	Random Init-R18	0.187 (0.154 - 0.226)	0.966 (0.958 - 0.973)	0.410 (0.355 - 0.463)	0.990 (0.989 - 0.990)	0.276 (0.242 - 0.312)	0.286 (0.249 - 0.322)
	ECGCLIP-R18	0.196 (0.165 - 0.237)	0.968 (0.959 - 0.975)	0.353 (0.304 - 0.405)	0.993 (0.992 - 0.994)	0.290 (0.251 - 0.333)	0.289 (0.249 - 0.333)
	ECGCLIP-R34	0.195 (0.161 - 0.237)	0.967 (0.958 - 0.975)	0.321 (0.272 - 0.373)	0.994 (0.993 - 0.995)	0.286 (0.244 - 0.327)	0.282 (0.241 - 0.324)
	Merl-R18	0.189 (0.157 - 0.228)	0.966 (0.957 - 0.974)	0.288 (0.237 - 0.341)	0.993 (0.993 - 0.994)	0.252 (0.208 - 0.296)	0.249 (0.204 - 0.293)
LVDD	Random Init-R18	0.035 (0.029 - 0.045)	0.776 (0.755 - 0.795)	0.181 (0.146 - 0.217)	0.967 (0.966 - 0.969)	0.078 (0.062 - 0.094)	0.079 (0.060 - 0.098)
	ECGCLIP-R18	0.043 (0.035 - 0.057)	0.793 (0.773 - 0.812)	0.129 (0.099 - 0.160)	0.986 (0.985 - 0.987)	0.099 (0.076 - 0.121)	0.091 (0.067 - 0.114)
	ECGCLIP-R34	0.045 (0.035 - 0.057)	0.799 (0.780 - 0.816)	0.255 (0.216 - 0.296)	0.958 (0.956 - 0.959)	0.090 (0.074 - 0.106)	0.100 (0.081 - 0.121)
	Merl-R18	0.039 (0.031 - 0.050)	0.785 (0.767 - 0.804)	0.181 (0.146 - 0.218)	0.972 (0.970 - 0.973)	0.087 (0.070 - 0.104)	0.086 (0.067 - 0.106)
VA	Random Init-R18	0.129 (0.097 - 0.160)	0.953 (0.935 - 0.967)	0.290 (0.218 - 0.356)	0.996 (0.995 - 0.996)	0.232 (0.173 - 0.283)	0.234 (0.174 - 0.286)
	ECGCLIP-R18	0.161 (0.123 - 0.205)	0.964 (0.949 - 0.976)	0.314 (0.246 - 0.382)	0.996 (0.995 - 0.996)	0.251 (0.197 - 0.302)	0.253 (0.199 - 0.305)
	ECGCLIP-R34	0.195 (0.146 - 0.249)	0.967 (0.953 - 0.979)	0.296 (0.232 - 0.359)	0.997 (0.996 - 0.997)	0.278 (0.219 - 0.335)	0.276 (0.217 - 0.333)
	Merl-R18	0.143 (0.109 - 0.184)	0.961 (0.945 - 0.974)	0.331 (0.262 - 0.402)	0.994 (0.994 - 0.995)	0.230 (0.181 - 0.277)	0.238 (0.187 - 0.287)
RAE	Random Init-R18	0.013 (0.009 - 0.020)	0.877 (0.849 - 0.902)	0.079 (0.030 - 0.139)	0.991 (0.990 - 0.992)	0.030 (0.011 - 0.051)	0.034 (0.010 - 0.062)
	ECGCLIP-R18	0.014 (0.010 - 0.019)	0.886 (0.857 - 0.909)	0.386 (0.293 - 0.478)	0.945 (0.943 - 0.947)	0.029 (0.020 - 0.038)	0.067 (0.047 - 0.087)
	ECGCLIP-R34	0.014 (0.010 - 0.022)	0.892 (0.868 - 0.912)	0.347 (0.252 - 0.442)	0.950 (0.948 - 0.952)	0.028 (0.019 - 0.038)	0.063 (0.042 - 0.083)
	Merl-R18	0.015 (0.010 - 0.022)	0.880 (0.850 - 0.906)	0.099 (0.040 - 0.160)	0.989 (0.988 - 0.990)	0.032 (0.013 - 0.052)	0.039 (0.013 - 0.065)
AM	Random Init-R18	0.014 (0.007 - 0.036)	0.752 (0.690 - 0.807)	0.074 (0.023 - 0.138)	0.996 (0.995 - 0.997)	0.044 (0.014 - 0.084)	0.046 (0.013 - 0.089)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAM	ECGCLIP-R18	0.029 (0.008 - 0.070)	0.755 (0.694 - 0.809)	0.025 (0.000 - 0.064)	1.000 (1.000 - 1.000)	0.041 (0.000 - 0.103)	0.053 (-0.001 - 0.136)
	ECGCLIP-R34	0.018 (0.009 - 0.040)	0.783 (0.723 - 0.834)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Merl-R18	0.018 (0.008 - 0.048)	0.748 (0.687 - 0.804)	0.086 (0.028 - 0.156)	0.998 (0.997 - 0.998)	0.069 (0.022 - 0.122)	0.068 (0.020 - 0.123)
	Random Init-R18	0.013 (0.003 - 0.067)	0.847 (0.788 - 0.902)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.002 (0.001 - 0.002)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.014 (0.005 - 0.038)	0.865 (0.797 - 0.923)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.002 (0.001 - 0.002)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.022 (0.008 - 0.075)	0.882 (0.809 - 0.943)	0.167 (0.049 - 0.292)	0.996 (0.995 - 0.996)	0.048 (0.016 - 0.086)	0.067 (0.020 - 0.117)
	Merl-R18	0.068 (0.007 - 0.152)	0.879 (0.813 - 0.934)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.002 (0.001 - 0.002)	0.000 (0.000 - 0.000)
	Random Init-R18	0.370 (0.359 - 0.381)	0.788 (0.783 - 0.793)	0.598 (0.586 - 0.610)	0.789 (0.785 - 0.793)	0.428 (0.419 - 0.437)	0.311 (0.302 - 0.322)
	ECGCLIP-R18	0.393 (0.382 - 0.405)	0.806 (0.801 - 0.811)	0.602 (0.591 - 0.614)	0.808 (0.804 - 0.812)	0.448 (0.440 - 0.457)	0.337 (0.327 - 0.347)
	ECGCLIP-R34	0.399 (0.388 - 0.411)	0.810 (0.805 - 0.814)	0.615 (0.604 - 0.626)	0.805 (0.801 - 0.809)	0.452 (0.443 - 0.461)	0.342 (0.332 - 0.353)
AVC	Merl-R18	0.383 (0.372 - 0.395)	0.798 (0.793 - 0.803)	0.629 (0.616 - 0.640)	0.784 (0.780 - 0.788)	0.441 (0.433 - 0.450)	0.329 (0.320 - 0.339)
	Random Init-R18	0.536 (0.521 - 0.553)	0.907 (0.902 - 0.911)	0.619 (0.605 - 0.633)	0.941 (0.939 - 0.943)	0.551 (0.539 - 0.563)	0.507 (0.495 - 0.520)
	ECGCLIP-R18	0.568 (0.552 - 0.585)	0.916 (0.911 - 0.920)	0.618 (0.604 - 0.633)	0.948 (0.946 - 0.950)	0.570 (0.558 - 0.581)	0.528 (0.515 - 0.540)
	ECGCLIP-R34	0.575 (0.560 - 0.592)	0.916 (0.912 - 0.920)	0.627 (0.613 - 0.642)	0.947 (0.945 - 0.949)	0.573 (0.562 - 0.585)	0.532 (0.519 - 0.544)
TR	Merl-R18	0.560 (0.545 - 0.576)	0.912 (0.908 - 0.917)	0.584 (0.569 - 0.600)	0.953 (0.951 - 0.955)	0.561 (0.548 - 0.573)	0.518 (0.505 - 0.531)
	Random Init-R18	0.426 (0.410 - 0.442)	0.899 (0.894 - 0.904)	0.612 (0.596 - 0.628)	0.926 (0.924 - 0.928)	0.477 (0.465 - 0.489)	0.440 (0.427 - 0.452)
	ECGCLIP-R18	0.461 (0.444 - 0.477)	0.910 (0.906 - 0.915)	0.650 (0.634 - 0.666)	0.926 (0.924 - 0.928)	0.500 (0.488 - 0.512)	0.467 (0.453 - 0.479)
MR							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AR	ECGCLIP-R34	0.463 (0.446 - 0.480)	0.913 (0.908 - 0.917)	0.627 (0.610 - 0.643)	0.933 (0.931 - 0.935)	0.505 (0.491 - 0.518)	0.469 (0.455 - 0.482)
	Merl-R18	0.450 (0.434 - 0.467)	0.906 (0.901 - 0.911)	0.607 (0.591 - 0.622)	0.933 (0.931 - 0.936)	0.493 (0.480 - 0.505)	0.455 (0.441 - 0.468)
	Random Init-R18	0.321 (0.299 - 0.343)	0.874 (0.866 - 0.882)	0.354 (0.332 - 0.376)	0.975 (0.974 - 0.977)	0.365 (0.345 - 0.384)	0.339 (0.318 - 0.359)
	ECGCLIP-R18	0.388 (0.364 - 0.412)	0.890 (0.882 - 0.897)	0.413 (0.389 - 0.435)	0.977 (0.975 - 0.978)	0.420 (0.399 - 0.441)	0.396 (0.374 - 0.417)
	ECGCLIP-R34	0.395 (0.371 - 0.418)	0.891 (0.884 - 0.898)	0.388 (0.365 - 0.409)	0.980 (0.979 - 0.981)	0.417 (0.395 - 0.437)	0.396 (0.373 - 0.416)
	Merl-R18	0.370 (0.346 - 0.393)	0.883 (0.875 - 0.891)	0.414 (0.391 - 0.438)	0.973 (0.972 - 0.975)	0.404 (0.384 - 0.425)	0.379 (0.358 - 0.400)
MVPLAC	Random Init-R18	0.081 (0.073 - 0.092)	0.803 (0.792 - 0.814)	0.180 (0.158 - 0.205)	0.964 (0.963 - 0.966)	0.133 (0.117 - 0.151)	0.111 (0.095 - 0.130)
	ECGCLIP-R18	0.088 (0.079 - 0.100)	0.815 (0.804 - 0.825)	0.256 (0.230 - 0.282)	0.951 (0.949 - 0.953)	0.153 (0.137 - 0.170)	0.138 (0.121 - 0.155)
	ECGCLIP-R34	0.089 (0.080 - 0.101)	0.820 (0.809 - 0.830)	0.216 (0.191 - 0.240)	0.959 (0.957 - 0.961)	0.146 (0.129 - 0.163)	0.127 (0.109 - 0.144)
	Merl-R18	0.083 (0.074 - 0.094)	0.806 (0.794 - 0.816)	0.223 (0.197 - 0.248)	0.957 (0.955 - 0.959)	0.145 (0.128 - 0.161)	0.126 (0.109 - 0.143)
MS	Random Init-R18	0.438 (0.402 - 0.482)	0.962 (0.956 - 0.968)	0.461 (0.425 - 0.497)	0.992 (0.991 - 0.993)	0.470 (0.437 - 0.502)	0.461 (0.429 - 0.495)
	ECGCLIP-R18	0.488 (0.449 - 0.532)	0.970 (0.964 - 0.975)	0.534 (0.497 - 0.570)	0.992 (0.991 - 0.992)	0.519 (0.490 - 0.551)	0.511 (0.482 - 0.543)
	ECGCLIP-R34	0.505 (0.468 - 0.547)	0.968 (0.962 - 0.973)	0.498 (0.463 - 0.534)	0.993 (0.993 - 0.994)	0.525 (0.493 - 0.557)	0.518 (0.487 - 0.552)
AS	Merl-R18	0.462 (0.424 - 0.503)	0.966 (0.961 - 0.972)	0.497 (0.459 - 0.530)	0.992 (0.991 - 0.993)	0.498 (0.467 - 0.529)	0.490 (0.459 - 0.521)
	Random Init-R18	0.220 (0.192 - 0.250)	0.926 (0.916 - 0.936)	0.364 (0.327 - 0.401)	0.986 (0.985 - 0.987)	0.306 (0.274 - 0.333)	0.298 (0.266 - 0.326)
	ECGCLIP-R18	0.293 (0.258 - 0.328)	0.936 (0.927 - 0.945)	0.432 (0.393 - 0.468)	0.985 (0.983 - 0.986)	0.341 (0.311 - 0.369)	0.337 (0.307 - 0.366)
	ECGCLIP-R34	0.302 (0.267 - 0.340)	0.939 (0.930 - 0.947)	0.331 (0.295 - 0.365)	0.992 (0.991 - 0.993)	0.348 (0.314 - 0.380)	0.340 (0.306 - 0.372)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
BAV	Merl-R18	0.248 (0.218 - 0.280)	0.930 (0.921 - 0.939)	0.417 (0.378 - 0.453)	0.984 (0.983 - 0.985)	0.325 (0.295 - 0.351)	0.322 (0.291 - 0.349)
	Random Init-R18	0.065 (0.053 - 0.082)	0.809 (0.789 - 0.829)	0.147 (0.113 - 0.177)	0.989 (0.988 - 0.990)	0.131 (0.102 - 0.158)	0.122 (0.093 - 0.150)
	ECGCLIP-R18	0.086 (0.068 - 0.109)	0.833 (0.814 - 0.851)	0.202 (0.166 - 0.240)	0.986 (0.985 - 0.987)	0.157 (0.130 - 0.183)	0.150 (0.123 - 0.177)
	ECGCLIP-R34	0.105 (0.083 - 0.132)	0.842 (0.822 - 0.860)	0.223 (0.187 - 0.261)	0.988 (0.987 - 0.989)	0.188 (0.159 - 0.216)	0.181 (0.151 - 0.210)
PS	Merl-R18	0.069 (0.056 - 0.086)	0.818 (0.798 - 0.837)	0.213 (0.177 - 0.251)	0.983 (0.982 - 0.984)	0.147 (0.122 - 0.172)	0.143 (0.116 - 0.169)
	Random Init-R18	0.075 (0.023 - 0.118)	0.911 (0.851 - 0.959)	0.082 (0.015 - 0.167)	0.999 (0.998 - 0.999)	0.068 (0.015 - 0.136)	0.068 (0.013 - 0.137)
	ECGCLIP-R18	0.066 (0.035 - 0.114)	0.924 (0.871 - 0.968)	0.347 (0.209 - 0.477)	0.996 (0.996 - 0.997)	0.145 (0.083 - 0.206)	0.176 (0.107 - 0.245)
	ECGCLIP-R34	0.095 (0.053 - 0.173)	0.933 (0.879 - 0.976)	0.102 (0.024 - 0.196)	0.999 (0.999 - 0.999)	0.108 (0.026 - 0.194)	0.107 (0.025 - 0.193)
EA	Merl-R18	0.010 (0.005 - 0.016)	0.889 (0.838 - 0.930)	0.184 (0.079 - 0.302)	0.988 (0.987 - 0.989)	0.029 (0.010 - 0.048)	0.050 (0.019 - 0.083)
	Random Init-R18	0.003 (0.001 - 0.005)	0.777 (0.660 - 0.878)	0.000 (0.000 - 0.000)	0.999 (0.999 - 1.000)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - 0.000)
	ECGCLIP-R18	0.105 (0.032 - 0.235)	0.920 (0.866 - 0.964)	0.111 (0.000 - 0.240)	1.000 (1.000 - 1.000)	0.162 (0.000 - 0.326)	0.182 (0.000 - 0.365)
	ECGCLIP-R34	0.253 (0.080 - 0.433)	0.902 (0.828 - 0.958)	0.370 (0.200 - 0.559)	0.999 (0.999 - 0.999)	0.253 (0.127 - 0.368)	0.266 (0.139 - 0.390)
PR	Merl-R18	0.004 (0.001 - 0.012)	0.704 (0.565 - 0.825)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Random Init-R18	0.051 (0.001 - 0.131)	0.751 (0.617 - 0.873)	0.095 (0.000 - 0.242)	0.998 (0.997 - 0.998)	0.031 (0.000 - 0.079)	0.041 (-0.001 - 0.100)
	ECGCLIP-R18	0.037 (0.002 - 0.126)	0.809 (0.693 - 0.900)	0.190 (0.042 - 0.381)	0.998 (0.997 - 0.998)	0.059 (0.014 - 0.117)	0.081 (0.018 - 0.158)
	ECGCLIP-R34	0.008 (0.002 - 0.021)	0.857 (0.770 - 0.930)	0.095 (0.000 - 0.250)	0.999 (0.998 - 0.999)	0.047 (0.000 - 0.118)	0.054 (-0.001 - 0.134)
	Merl-R18	0.057 (0.002 - 0.151)	0.871 (0.814 - 0.928)	0.238 (0.056 - 0.445)	0.993 (0.992 - 0.994)	0.029 (0.006 - 0.057)	0.059 (0.012 - 0.112)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
TS	Random Init-R18	0.005 (0.001 - 0.029)	0.891 (0.807 - 0.960)	0.583 (0.285 - 0.857)	0.924 (0.922 - 0.927)	0.004 (0.001 - 0.007)	0.031 (0.011 - 0.049)
	ECGCLIP-R18	0.004 (0.001 - 0.012)	0.722 (0.507 - 0.917)	0.167 (0.000 - 0.429)	0.995 (0.994 - 0.995)	0.015 (0.000 - 0.037)	0.035 (-0.001 - 0.084)
	ECGCLIP-R34	0.003 (0.001 - 0.007)	0.856 (0.745 - 0.959)	0.083 (0.000 - 0.286)	0.995 (0.994 - 0.995)	0.007 (0.000 - 0.028)	0.017 (-0.001 - 0.061)
	Merl-R18	0.001 (0.000 - 0.001)	0.681 (0.489 - 0.854)	0.000 (0.000 - 0.000)	0.999 (0.999 - 0.999)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
CHD	Random Init-R18	0.343 (0.308 - 0.379)	0.890 (0.876 - 0.903)	0.427 (0.392 - 0.460)	0.990 (0.989 - 0.990)	0.427 (0.395 - 0.455)	0.416 (0.385 - 0.445)
	ECGCLIP-R18	0.384 (0.347 - 0.420)	0.901 (0.887 - 0.914)	0.494 (0.458 - 0.529)	0.989 (0.988 - 0.990)	0.471 (0.440 - 0.500)	0.461 (0.430 - 0.491)
	ECGCLIP-R34	0.411 (0.372 - 0.444)	0.908 (0.895 - 0.920)	0.490 (0.454 - 0.525)	0.990 (0.989 - 0.990)	0.474 (0.443 - 0.501)	0.464 (0.433 - 0.493)
	Merl-R18	0.361 (0.324 - 0.398)	0.897 (0.883 - 0.910)	0.480 (0.445 - 0.513)	0.988 (0.987 - 0.989)	0.454 (0.424 - 0.482)	0.444 (0.414 - 0.472)
ASD	Random Init-R18	0.526 (0.494 - 0.560)	0.920 (0.908 - 0.931)	0.477 (0.446 - 0.510)	0.993 (0.993 - 0.994)	0.528 (0.500 - 0.558)	0.523 (0.495 - 0.554)
	ECGCLIP-R18	0.599 (0.566 - 0.631)	0.931 (0.920 - 0.941)	0.549 (0.517 - 0.582)	0.993 (0.993 - 0.994)	0.585 (0.557 - 0.613)	0.578 (0.551 - 0.607)
	ECGCLIP-R34	0.619 (0.588 - 0.650)	0.938 (0.928 - 0.947)	0.606 (0.576 - 0.637)	0.991 (0.991 - 0.992)	0.598 (0.573 - 0.624)	0.589 (0.564 - 0.616)
	Merl-R18	0.566 (0.534 - 0.601)	0.928 (0.917 - 0.938)	0.568 (0.536 - 0.600)	0.990 (0.989 - 0.991)	0.556 (0.530 - 0.584)	0.547 (0.520 - 0.575)
VSD	Random Init-R18	0.174 (0.132 - 0.218)	0.868 (0.847 - 0.891)	0.178 (0.139 - 0.219)	0.996 (0.995 - 0.996)	0.203 (0.160 - 0.245)	0.200 (0.157 - 0.242)
	ECGCLIP-R18	0.203 (0.160 - 0.250)	0.889 (0.869 - 0.908)	0.317 (0.268 - 0.367)	0.992 (0.992 - 0.993)	0.274 (0.232 - 0.317)	0.270 (0.228 - 0.314)
	ECGCLIP-R34	0.236 (0.191 - 0.284)	0.903 (0.885 - 0.920)	0.312 (0.264 - 0.362)	0.995 (0.995 - 0.996)	0.318 (0.273 - 0.363)	0.313 (0.268 - 0.358)
	Merl-R18	0.176 (0.134 - 0.222)	0.873 (0.850 - 0.895)	0.289 (0.242 - 0.338)	0.991 (0.990 - 0.992)	0.235 (0.196 - 0.273)	0.232 (0.194 - 0.272)
PFO	Random Init-R18	0.016 (0.012 - 0.022)	0.766 (0.730 - 0.800)	0.162 (0.111 - 0.212)	0.975 (0.973 - 0.976)	0.045 (0.031 - 0.061)	0.056 (0.035 - 0.077)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
PDA	ECGCLIP-R18	0.018 (0.014 - 0.025)	0.780 (0.745 - 0.815)	0.192 (0.140 - 0.251)	0.973 (0.971 - 0.974)	0.051 (0.035 - 0.066)	0.065 (0.043 - 0.089)
	ECGCLIP-R34	0.022 (0.016 - 0.030)	0.791 (0.756 - 0.826)	0.268 (0.203 - 0.332)	0.966 (0.964 - 0.967)	0.057 (0.043 - 0.073)	0.082 (0.059 - 0.105)
	Merl-R18	0.019 (0.014 - 0.027)	0.785 (0.752 - 0.819)	0.177 (0.124 - 0.233)	0.973 (0.972 - 0.974)	0.047 (0.032 - 0.062)	0.059 (0.038 - 0.081)
	Random Init-R18	0.079 (0.046 - 0.135)	0.860 (0.824 - 0.891)	0.187 (0.118 - 0.255)	0.998 (0.997 - 0.998)	0.176 (0.112 - 0.235)	0.174 (0.110 - 0.235)
	ECGCLIP-R18	0.111 (0.062 - 0.173)	0.881 (0.849 - 0.911)	0.187 (0.121 - 0.257)	0.998 (0.998 - 0.999)	0.202 (0.133 - 0.272)	0.200 (0.131 - 0.271)
	ECGCLIP-R34	0.113 (0.064 - 0.175)	0.893 (0.862 - 0.921)	0.146 (0.087 - 0.208)	0.999 (0.998 - 0.999)	0.181 (0.108 - 0.250)	0.185 (0.111 - 0.256)
	Merl-R18	0.092 (0.049 - 0.151)	0.866 (0.834 - 0.896)	0.195 (0.129 - 0.271)	0.997 (0.996 - 0.997)	0.166 (0.108 - 0.225)	0.165 (0.108 - 0.225)
	Random Init-R18	0.000 (0.000 - 0.002)	0.633 (0.275 - 0.988)	0.500 (0.000 - 1.000)	0.952 (0.950 - 0.954)	0.001 (0.000 - 0.003)	0.014 (-0.002 - 0.035)
	ECGCLIP-R18	0.001 (0.000 - 0.002)	0.970 (0.964 - 0.976)	1.000 (1.000 - 1.000)	0.872 (0.869 - 0.875)	0.001 (0.000 - 0.002)	0.017 (0.012 - 0.027)
	ECGCLIP-R34	0.001 (0.000 - 0.003)	0.980 (0.974 - 0.987)	1.000 (1.000 - 1.000)	0.865 (0.862 - 0.868)	0.001 (0.000 - 0.002)	0.016 (0.012 - 0.026)
TOF	Merl-R18	0.000 (0.000 - 0.000)	0.413 (0.073 - 0.755)	0.000 (0.000 - 0.000)	0.794 (0.791 - 0.798)	0.000 (0.000 - 0.000)	-0.003 (-0.005 - -0.002)
	Random Init-R18	0.473 (0.461 - 0.486)	0.842 (0.837 - 0.846)	0.649 (0.638 - 0.660)	0.850 (0.847 - 0.854)	0.509 (0.500 - 0.518)	0.421 (0.411 - 0.431)
	ECGCLIP-R18	0.498 (0.485 - 0.511)	0.852 (0.847 - 0.857)	0.606 (0.593 - 0.617)	0.882 (0.878 - 0.885)	0.523 (0.513 - 0.532)	0.436 (0.426 - 0.446)
	ECGCLIP-R34	0.502 (0.489 - 0.515)	0.853 (0.848 - 0.858)	0.596 (0.584 - 0.608)	0.887 (0.884 - 0.890)	0.523 (0.514 - 0.532)	0.437 (0.426 - 0.448)
PH	Merl-R18	0.486 (0.473 - 0.499)	0.846 (0.841 - 0.851)	0.637 (0.625 - 0.648)	0.858 (0.855 - 0.862)	0.511 (0.502 - 0.521)	0.423 (0.412 - 0.434)
	Random Init-R18	0.273 (0.261 - 0.287)	0.812 (0.805 - 0.819)	0.443 (0.427 - 0.459)	0.902 (0.899 - 0.905)	0.349 (0.337 - 0.361)	0.285 (0.272 - 0.298)
	ECGCLIP-R18	0.302 (0.288 - 0.316)	0.830 (0.824 - 0.836)	0.491 (0.475 - 0.506)	0.898 (0.895 - 0.901)	0.373 (0.361 - 0.385)	0.313 (0.301 - 0.326)
AD							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
Dex	ECGCLIP-R34	0.305 (0.291 - 0.319)	0.833 (0.826 - 0.839)	0.551 (0.535 - 0.566)	0.877 (0.873 - 0.879)	0.376 (0.364 - 0.387)	0.321 (0.309 - 0.334)
	Merl-R18	0.287 (0.273 - 0.301)	0.822 (0.816 - 0.828)	0.508 (0.492 - 0.525)	0.886 (0.883 - 0.889)	0.365 (0.353 - 0.377)	0.306 (0.294 - 0.319)
	Random Init-R18	0.000 (0.000 - 0.001)	0.622 (0.486 - 0.754)	0.000 (0.000 - 0.000)	0.986 (0.985 - 0.987)	0.000 (0.000 - 0.000)	-0.002 (-0.003 - -0.002)
	ECGCLIP-R18	0.084 (0.002 - 0.235)	0.853 (0.697 - 0.973)	0.333 (0.091 - 0.589)	0.990 (0.989 - 0.991)	0.020 (0.004 - 0.039)	0.057 (0.013 - 0.102)
	ECGCLIP-R34	0.121 (0.004 - 0.312)	0.832 (0.676 - 0.968)	0.200 (0.000 - 0.455)	0.999 (0.999 - 0.999)	0.085 (0.000 - 0.184)	0.103 (-0.001 - 0.216)
	Merl-R18	0.004 (0.001 - 0.014)	0.746 (0.572 - 0.909)	0.067 (0.000 - 0.222)	0.998 (0.998 - 0.998)	0.019 (0.000 - 0.060)	0.026 (-0.001 - 0.084)
PE	Random Init-R18	0.280 (0.264 - 0.300)	0.852 (0.845 - 0.860)	0.424 (0.406 - 0.444)	0.946 (0.944 - 0.948)	0.349 (0.334 - 0.365)	0.314 (0.297 - 0.330)
	ECGCLIP-R18	0.310 (0.293 - 0.330)	0.866 (0.859 - 0.873)	0.489 (0.470 - 0.510)	0.941 (0.939 - 0.943)	0.377 (0.363 - 0.393)	0.346 (0.331 - 0.363)
	ECGCLIP-R34	0.318 (0.300 - 0.339)	0.867 (0.860 - 0.874)	0.458 (0.438 - 0.478)	0.950 (0.948 - 0.952)	0.384 (0.368 - 0.401)	0.351 (0.334 - 0.368)
	Merl-R18	0.293 (0.277 - 0.313)	0.859 (0.852 - 0.866)	0.475 (0.455 - 0.495)	0.938 (0.936 - 0.940)	0.361 (0.346 - 0.377)	0.329 (0.314 - 0.345)
	Random Init-R18	0.043 (0.027 - 0.070)	0.959 (0.940 - 0.974)	0.281 (0.169 - 0.403)	0.994 (0.993 - 0.995)	0.092 (0.051 - 0.137)	0.122 (0.071 - 0.179)
CP	ECGCLIP-R18	0.111 (0.058 - 0.198)	0.976 (0.965 - 0.986)	0.246 (0.137 - 0.362)	0.999 (0.998 - 0.999)	0.217 (0.127 - 0.308)	0.217 (0.127 - 0.310)
	ECGCLIP-R34	0.175 (0.096 - 0.283)	0.979 (0.965 - 0.990)	0.386 (0.267 - 0.510)	0.998 (0.998 - 0.999)	0.278 (0.189 - 0.370)	0.289 (0.200 - 0.380)
	Merl-R18	0.037 (0.023 - 0.063)	0.964 (0.951 - 0.975)	0.140 (0.060 - 0.232)	0.997 (0.997 - 0.998)	0.081 (0.033 - 0.135)	0.088 (0.036 - 0.147)
AF	Random Init-R18	0.953 (0.949 - 0.956)	0.994 (0.994 - 0.995)	0.943 (0.939 - 0.947)	0.985 (0.984 - 0.986)	0.921 (0.917 - 0.924)	0.909 (0.905 - 0.913)
	ECGCLIP-R18	0.965 (0.961 - 0.968)	0.996 (0.995 - 0.996)	0.951 (0.947 - 0.954)	0.987 (0.986 - 0.988)	0.932 (0.929 - 0.935)	0.922 (0.919 - 0.926)
	ECGCLIP-R34	0.970 (0.967 - 0.974)	0.996 (0.996 - 0.997)	0.956 (0.952 - 0.959)	0.988 (0.987 - 0.989)	0.938 (0.935 - 0.941)	0.929 (0.926 - 0.932)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
	Merl-R18	0.967 (0.964 - 0.971)	0.996 (0.995 - 0.996)	0.950 (0.946 - 0.953)	0.987 (0.986 - 0.988)	0.931 (0.928 - 0.934)	0.921 (0.918 - 0.925)
Tier 3: Rare Diseases							
CA	Random Init-R18	0.077 (0.038 - 0.087)	0.974 (0.962 - 0.983)	0.061 (0.019 - 0.106)	0.999 (0.999 - 0.999)	0.065 (0.020 - 0.113)	0.064 (0.019 - 0.112)
	ECGCLIP-R18	0.165 (0.109 - 0.232)	0.971 (0.957 - 0.981)	0.270 (0.192 - 0.350)	0.999 (0.999 - 0.999)	0.252 (0.182 - 0.327)	0.252 (0.182 - 0.327)
	ECGCLIP-R34	0.201 (0.143 - 0.275)	0.981 (0.970 - 0.990)	0.252 (0.179 - 0.333)	0.999 (0.999 - 1.000)	0.276 (0.200 - 0.356)	0.277 (0.200 - 0.357)
ARVC	Merl-R18	0.049 (0.028 - 0.081)	0.923 (0.888 - 0.952)	0.139 (0.077 - 0.208)	0.999 (0.998 - 0.999)	0.118 (0.064 - 0.173)	0.118 (0.064 - 0.176)
	Random Init-R18	0.000 (0.000 - 0.001)	0.940 (0.856 - 0.987)	0.000 (0.000 - 0.000)	0.995 (0.995 - 0.996)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
	ECGCLIP-R18	0.006 (0.001 - 0.018)	0.997 (0.995 - 0.999)	0.250 (0.000 - 1.000)	0.999 (0.999 - 1.000)	0.028 (0.000 - 0.085)	0.061 (0.000 - 0.166)
	ECGCLIP-R34	0.019 (0.004 - 0.056)	0.999 (0.999 - 1.000)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Merl-R18	0.004 (0.000 - 0.023)	0.990 (0.981 - 0.999)	0.250 (0.000 - 0.750)	0.999 (0.998 - 0.999)	0.013 (0.000 - 0.045)	0.041 (0.000 - 0.122)
TC	Random Init-R18	0.000 (0.000 - 0.001)	0.789 (0.709 - 0.870)	0.929 (0.750 - 1.000)	0.519 (0.516 - 0.522)	0.000 (0.000 - 0.001)	0.010 (0.005 - 0.014)
	ECGCLIP-R18	0.001 (0.000 - 0.003)	0.814 (0.648 - 0.925)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.002 (0.001 - 0.006)	0.898 (0.744 - 0.973)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Merl-R18	0.000 (0.000 - 0.001)	0.735 (0.602 - 0.860)	0.000 (0.000 - 0.000)	0.993 (0.993 - 0.994)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)
AC	Random Init-R18	0.001 (0.000 - 0.007)	0.695 (0.465 - 0.902)	0.167 (0.000 - 0.517)	0.959 (0.958 - 0.960)	0.000 (0.000 - 0.001)	0.005 (-0.002 - 0.017)
	ECGCLIP-R18	0.001 (0.000 - 0.002)	0.935 (0.848 - 0.986)	0.000 (0.000 - 0.000)	0.999 (0.998 - 0.999)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.001 (0.000 - 0.001)	0.882 (0.749 - 0.980)	0.000 (0.000 - 0.000)	0.991 (0.990 - 0.991)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - 0.000)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
NCVM	Merl-R18	0.001 (0.000 - 0.003)	0.762 (0.477 - 0.970)	0.000 (0.000 - 0.000)	0.998 (0.998 - 0.998)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Random Init-R18	0.000 (0.000 - 0.000)	0.470 (0.327 - 0.622)	0.438 (0.182 - 0.688)	0.494 (0.492 - 0.497)	0.000 (0.000 - 0.000)	-0.002 (-0.007 - 0.004)
	ECGCLIP-R18	0.000 (0.000 - 0.001)	0.765 (0.678 - 0.839)	0.063 (0.000 - 0.200)	0.982 (0.981 - 0.983)	0.001 (0.000 - 0.003)	0.004 (-0.002 - 0.017)
	ECGCLIP-R34	0.000 (0.000 - 0.000)	0.572 (0.480 - 0.674)	0.000 (0.000 - 0.000)	0.984 (0.984 - 0.985)	0.000 (0.000 - 0.000)	-0.002 (-0.002 - -0.001)
	Merl-R18	0.000 (0.000 - 0.001)	0.515 (0.339 - 0.680)	0.250 (0.055 - 0.500)	0.804 (0.801 - 0.806)	0.000 (0.000 - 0.001)	0.002 (-0.004 - 0.008)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S19: Detailed task-specific performance comparison on the external Xiamen cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
OMI	Random Init-R18	0.172 (0.116 – 0.266)	0.994 (0.990 – 0.997)	0.702 (0.558 – 0.833)	0.995 (0.995 – 0.996)	0.227 (0.162 – 0.294)	0.307 (0.234 – 0.374)
	ECGCLIP-R18	0.196 (0.130 – 0.299)	0.997 (0.995 – 0.998)	0.851 (0.750 – 0.953)	0.994 (0.993 – 0.994)	0.220 (0.164 – 0.280)	0.327 (0.267 – 0.386)
	ECGCLIP-R34	0.224 (0.145 – 0.336)	0.997 (0.996 – 0.998)	0.809 (0.696 – 0.915)	0.994 (0.993 – 0.994)	0.210 (0.154 – 0.268)	0.311 (0.251 – 0.369)
	Merl-R18	0.208 (0.139 – 0.302)	0.997 (0.995 – 0.998)	0.851 (0.743 – 0.950)	0.994 (0.993 – 0.994)	0.222 (0.166 – 0.280)	0.328 (0.269 – 0.386)
STEMI	Random Init-R18	0.031 (0.020 – 0.047)	0.717 (0.698 – 0.734)	0.028 (0.014 – 0.045)	0.998 (0.998 – 0.999)	0.047 (0.024 – 0.075)	0.061 (0.030 – 0.099)
	ECGCLIP-R18	0.039 (0.027 – 0.059)	0.756 (0.740 – 0.773)	0.035 (0.019 – 0.053)	0.998 (0.998 – 0.999)	0.057 (0.031 – 0.085)	0.071 (0.038 – 0.107)
	ECGCLIP-R34	0.064 (0.052 – 0.082)	0.893 (0.883 – 0.902)	0.044 (0.026 – 0.065)	0.998 (0.998 – 0.998)	0.071 (0.042 – 0.104)	0.084 (0.049 – 0.124)
	Merl-R18	0.024 (0.015 – 0.040)	0.638 (0.617 – 0.659)	0.032 (0.016 – 0.050)	0.998 (0.997 – 0.998)	0.052 (0.027 – 0.079)	0.059 (0.029 – 0.093)
NSR	Random Init-R18	0.967 (0.964 – 0.969)	0.926 (0.922 – 0.929)	0.984 (0.982 – 0.985)	0.725 (0.716 – 0.733)	0.953 (0.952 – 0.955)	0.777 (0.769 – 0.784)
	ECGCLIP-R18	0.970 (0.968 – 0.972)	0.932 (0.929 – 0.936)	0.980 (0.979 – 0.982)	0.759 (0.751 – 0.767)	0.956 (0.954 – 0.957)	0.793 (0.786 – 0.799)
	ECGCLIP-R34	0.977 (0.976 – 0.979)	0.945 (0.942 – 0.947)	0.989 (0.988 – 0.990)	0.749 (0.740 – 0.757)	0.959 (0.958 – 0.961)	0.807 (0.801 – 0.814)
	Merl-R18	0.967 (0.964 – 0.969)	0.927 (0.923 – 0.930)	0.979 (0.978 – 0.981)	0.744 (0.735 – 0.752)	0.954 (0.952 – 0.955)	0.781 (0.773 – 0.787)
SBrad	Random Init-R18	0.961 (0.956 – 0.965)	0.995 (0.994 – 0.995)	0.910 (0.903 – 0.918)	0.987 (0.985 – 0.988)	0.907 (0.901 – 0.913)	0.894 (0.888 – 0.900)
	ECGCLIP-R18	0.967 (0.962 – 0.970)	0.995 (0.995 – 0.996)	0.919 (0.912 – 0.926)	0.988 (0.987 – 0.989)	0.917 (0.912 – 0.922)	0.905 (0.899 – 0.911)
	ECGCLIP-R34	0.976 (0.973 – 0.979)	0.997 (0.996 – 0.997)	0.934 (0.928 – 0.941)	0.991 (0.990 – 0.992)	0.935 (0.931 – 0.940)	0.926 (0.921 – 0.932)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STach	Merl-R18	0.958 (0.953 - 0.963)	0.994 (0.993 - 0.995)	0.898 (0.889 - 0.905)	0.987 (0.986 - 0.988)	0.901 (0.895 - 0.907)	0.888 (0.881 - 0.895)
	Random Init-R18	0.964 (0.952 - 0.974)	0.999 (0.999 - 0.999)	0.918 (0.902 - 0.931)	0.998 (0.998 - 0.998)	0.926 (0.917 - 0.936)	0.924 (0.914 - 0.934)
	ECGCLIP-R18	0.969 (0.960 - 0.977)	0.999 (0.999 - 0.999)	0.938 (0.925 - 0.951)	0.998 (0.998 - 0.998)	0.938 (0.929 - 0.947)	0.936 (0.927 - 0.946)
	ECGCLIP-R34	0.974 (0.965 - 0.982)	0.999 (0.999 - 0.999)	0.954 (0.942 - 0.964)	0.998 (0.998 - 0.999)	0.948 (0.939 - 0.957)	0.946 (0.937 - 0.955)
VPB	Merl-R18	0.964 (0.954 - 0.973)	0.999 (0.999 - 0.999)	0.929 (0.914 - 0.941)	0.997 (0.997 - 0.998)	0.925 (0.915 - 0.935)	0.923 (0.912 - 0.933)
	Random Init-R18	0.672 (0.643 - 0.698)	0.975 (0.971 - 0.979)	0.695 (0.671 - 0.720)	0.993 (0.992 - 0.994)	0.719 (0.700 - 0.738)	0.712 (0.692 - 0.731)
	ECGCLIP-R18	0.715 (0.687 - 0.741)	0.983 (0.980 - 0.986)	0.768 (0.744 - 0.791)	0.993 (0.992 - 0.994)	0.766 (0.748 - 0.783)	0.759 (0.741 - 0.776)
	ECGCLIP-R34	0.767 (0.740 - 0.790)	0.992 (0.990 - 0.993)	0.872 (0.852 - 0.889)	0.992 (0.992 - 0.993)	0.819 (0.803 - 0.834)	0.815 (0.799 - 0.830)
APB	Merl-R18	0.689 (0.662 - 0.716)	0.979 (0.975 - 0.983)	0.716 (0.692 - 0.741)	0.993 (0.992 - 0.994)	0.734 (0.715 - 0.751)	0.726 (0.707 - 0.744)
	Random Init-R18	0.135 (0.122 - 0.152)	0.850 (0.837 - 0.860)	0.345 (0.315 - 0.373)	0.960 (0.958 - 0.962)	0.226 (0.205 - 0.244)	0.216 (0.194 - 0.234)
	ECGCLIP-R18	0.181 (0.162 - 0.203)	0.875 (0.863 - 0.884)	0.342 (0.312 - 0.370)	0.974 (0.973 - 0.976)	0.280 (0.256 - 0.302)	0.265 (0.240 - 0.287)
	ECGCLIP-R34	0.360 (0.329 - 0.391)	0.928 (0.919 - 0.936)	0.487 (0.454 - 0.518)	0.985 (0.984 - 0.986)	0.460 (0.433 - 0.484)	0.447 (0.420 - 0.472)
SArr	Merl-R18	0.151 (0.134 - 0.169)	0.858 (0.847 - 0.868)	0.336 (0.307 - 0.363)	0.969 (0.967 - 0.970)	0.251 (0.229 - 0.272)	0.237 (0.216 - 0.258)
	Random Init-R18	0.097 (0.088 - 0.108)	0.747 (0.736 - 0.759)	0.313 (0.289 - 0.337)	0.916 (0.913 - 0.918)	0.170 (0.156 - 0.184)	0.144 (0.129 - 0.159)
	ECGCLIP-R18	0.108 (0.098 - 0.119)	0.759 (0.747 - 0.771)	0.424 (0.400 - 0.452)	0.877 (0.874 - 0.880)	0.174 (0.163 - 0.186)	0.162 (0.149 - 0.176)
	ECGCLIP-R34	0.113 (0.103 - 0.125)	0.771 (0.760 - 0.782)	0.438 (0.414 - 0.465)	0.877 (0.874 - 0.880)	0.179 (0.167 - 0.191)	0.169 (0.155 - 0.183)
	Merl-R18	0.098 (0.089 - 0.108)	0.750 (0.738 - 0.763)	0.298 (0.275 - 0.321)	0.922 (0.919 - 0.924)	0.170 (0.156 - 0.184)	0.143 (0.127 - 0.158)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AF	Random Init-R18	0.784 (0.743 - 0.835)	0.998 (0.997 - 0.998)	0.906 (0.879 - 0.931)	0.996 (0.996 - 0.997)	0.794 (0.768 - 0.821)	0.798 (0.773 - 0.824)
	ECGCLIP-R18	0.821 (0.783 - 0.862)	0.998 (0.997 - 0.999)	0.895 (0.864 - 0.923)	0.997 (0.996 - 0.997)	0.803 (0.778 - 0.830)	0.805 (0.780 - 0.831)
	ECGCLIP-R34	0.874 (0.837 - 0.910)	0.998 (0.997 - 0.999)	0.948 (0.927 - 0.967)	0.997 (0.997 - 0.998)	0.848 (0.825 - 0.871)	0.851 (0.829 - 0.874)
	Merl-R18	0.797 (0.756 - 0.845)	0.998 (0.997 - 0.999)	0.909 (0.879 - 0.935)	0.997 (0.996 - 0.997)	0.802 (0.777 - 0.830)	0.806 (0.780 - 0.832)
AFL	Random Init-R18	0.642 (0.548 - 0.742)	0.997 (0.996 - 0.999)	0.614 (0.520 - 0.709)	0.999 (0.999 - 0.999)	0.614 (0.533 - 0.689)	0.613 (0.533 - 0.688)
	ECGCLIP-R18	0.700 (0.612 - 0.788)	0.999 (0.998 - 0.999)	0.675 (0.585 - 0.755)	0.999 (0.999 - 0.999)	0.653 (0.578 - 0.719)	0.652 (0.579 - 0.719)
	ECGCLIP-R34	0.725 (0.644 - 0.801)	0.998 (0.997 - 0.999)	0.711 (0.627 - 0.789)	0.999 (0.999 - 0.999)	0.701 (0.629 - 0.763)	0.701 (0.628 - 0.763)
	Merl-R18	0.649 (0.554 - 0.745)	0.998 (0.997 - 0.999)	0.658 (0.566 - 0.748)	0.999 (0.999 - 0.999)	0.636 (0.560 - 0.700)	0.635 (0.560 - 0.700)
JPB	Random Init-R18	0.023 (0.019 - 0.029)	0.819 (0.796 - 0.843)	0.096 (0.063 - 0.135)	0.984 (0.982 - 0.985)	0.048 (0.031 - 0.067)	0.046 (0.027 - 0.068)
	ECGCLIP-R18	0.031 (0.024 - 0.039)	0.853 (0.831 - 0.874)	0.060 (0.032 - 0.089)	0.991 (0.991 - 0.992)	0.047 (0.025 - 0.069)	0.041 (0.019 - 0.064)
	ECGCLIP-R34	0.080 (0.063 - 0.102)	0.920 (0.904 - 0.935)	0.164 (0.119 - 0.211)	0.993 (0.992 - 0.994)	0.139 (0.102 - 0.176)	0.135 (0.098 - 0.173)
	Merl-R18	0.025 (0.020 - 0.032)	0.834 (0.811 - 0.857)	0.040 (0.018 - 0.065)	0.994 (0.993 - 0.995)	0.039 (0.016 - 0.063)	0.033 (0.011 - 0.058)
ATach	Random Init-R18	0.080 (0.057 - 0.111)	0.944 (0.930 - 0.958)	0.248 (0.182 - 0.322)	0.995 (0.994 - 0.996)	0.177 (0.126 - 0.229)	0.181 (0.131 - 0.235)
	ECGCLIP-R18	0.129 (0.090 - 0.180)	0.951 (0.935 - 0.965)	0.262 (0.188 - 0.331)	0.996 (0.996 - 0.997)	0.223 (0.163 - 0.278)	0.222 (0.162 - 0.278)
	ECGCLIP-R34	0.153 (0.112 - 0.208)	0.959 (0.944 - 0.972)	0.338 (0.257 - 0.416)	0.996 (0.995 - 0.996)	0.258 (0.198 - 0.317)	0.262 (0.201 - 0.322)
	Merl-R18	0.090 (0.062 - 0.126)	0.941 (0.924 - 0.956)	0.214 (0.145 - 0.279)	0.996 (0.996 - 0.997)	0.187 (0.128 - 0.239)	0.185 (0.125 - 0.239)
JTach	Random Init-R18	0.273 (0.191 - 0.366)	0.973 (0.962 - 0.984)	0.158 (0.092 - 0.227)	1.000 (0.999 - 1.000)	0.234 (0.143 - 0.323)	0.265 (0.167 - 0.359)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
JEB	ECGCLIP-R18	0.360 (0.268 - 0.454)	0.979 (0.967 - 0.987)	0.246 (0.170 - 0.327)	1.000 (1.000 - 1.000)	0.361 (0.259 - 0.453)	0.409 (0.309 - 0.498)
	ECGCLIP-R34	0.382 (0.295 - 0.475)	0.981 (0.973 - 0.988)	0.211 (0.141 - 0.290)	1.000 (1.000 - 1.000)	0.327 (0.230 - 0.426)	0.391 (0.292 - 0.482)
	Merl-R18	0.285 (0.202 - 0.374)	0.972 (0.959 - 0.982)	0.211 (0.140 - 0.286)	1.000 (0.999 - 1.000)	0.312 (0.217 - 0.400)	0.355 (0.258 - 0.444)
	Random Init-R18	0.038 (0.008 - 0.184)	0.982 (0.971 - 0.991)	0.125 (0.000 - 0.300)	0.999 (0.999 - 0.999)	0.066 (0.000 - 0.154)	0.074 (-0.001 - 0.176)
	ECGCLIP-R18	0.081 (0.007 - 0.231)	0.979 (0.967 - 0.990)	0.125 (0.000 - 0.300)	0.999 (0.999 - 1.000)	0.089 (0.000 - 0.207)	0.092 (-0.001 - 0.218)
	ECGCLIP-R34	0.083 (0.009 - 0.221)	0.987 (0.980 - 0.993)	0.063 (0.000 - 0.200)	1.000 (0.999 - 1.000)	0.061 (0.000 - 0.182)	0.060 (0.000 - 0.184)
	Merl-R18	0.050 (0.007 - 0.227)	0.977 (0.963 - 0.990)	0.125 (0.000 - 0.300)	0.999 (0.999 - 1.000)	0.095 (0.000 - 0.216)	0.098 (0.000 - 0.222)
	Random Init-R18	0.797 (0.595 - 0.948)	1.000 (0.999 - 1.000)	0.688 (0.454 - 0.909)	1.000 (1.000 - 1.000)	0.759 (0.556 - 0.909)	0.763 (0.570 - 0.913)
	ECGCLIP-R18	0.840 (0.665 - 0.974)	1.000 (0.999 - 1.000)	0.750 (0.526 - 0.941)	1.000 (1.000 - 1.000)	0.828 (0.666 - 0.957)	0.832 (0.667 - 0.957)
	ECGCLIP-R34	0.800 (0.597 - 0.957)	0.999 (0.998 - 1.000)	0.688 (0.454 - 0.909)	1.000 (1.000 - 1.000)	0.786 (0.588 - 0.929)	0.794 (0.612 - 0.931)
SVT	Merl-R18	0.791 (0.582 - 0.947)	0.999 (0.998 - 1.000)	0.750 (0.526 - 0.941)	1.000 (1.000 - 1.000)	0.828 (0.640 - 0.957)	0.832 (0.666 - 0.957)
	Random Init-R18	0.024 (0.009 - 0.050)	0.968 (0.915 - 0.995)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.038 (0.012 - 0.092)	0.971 (0.923 - 0.995)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.058 (0.017 - 0.152)	0.979 (0.952 - 0.996)	0.067 (0.000 - 0.222)	1.000 (1.000 - 1.000)	0.100 (0.000 - 0.300)	0.115 (0.000 - 0.343)
VEB	Merl-R18	0.025 (0.009 - 0.054)	0.964 (0.907 - 0.994)	0.067 (0.000 - 0.200)	0.999 (0.999 - 1.000)	0.045 (0.000 - 0.133)	0.047 (-0.001 - 0.144)
	Random Init-R18	0.017 (0.004 - 0.047)	0.985 (0.972 - 0.994)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.025 (0.005 - 0.084)	0.991 (0.984 - 0.996)	0.091 (0.000 - 0.300)	1.000 (1.000 - 1.000)	0.083 (0.000 - 0.250)	0.083 (0.000 - 0.257)
VTach							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RBBB	ECGCLIP-R34	0.032 (0.006 - 0.185)	0.992 (0.988 - 0.995)	0.091 (0.000 - 0.300)	1.000 (1.000 - 1.000)	0.091 (0.000 - 0.296)	0.091 (0.000 - 0.300)
	Merl-R18	0.015 (0.004 - 0.044)	0.987 (0.980 - 0.994)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Random Init-R18	0.971 (0.963 - 0.978)	0.999 (0.999 - 0.999)	0.922 (0.907 - 0.935)	0.998 (0.997 - 0.998)	0.922 (0.911 - 0.933)	0.919 (0.908 - 0.931)
	ECGCLIP-R18	0.974 (0.967 - 0.981)	0.999 (0.999 - 0.999)	0.922 (0.907 - 0.936)	0.998 (0.998 - 0.998)	0.926 (0.916 - 0.937)	0.924 (0.913 - 0.935)
	ECGCLIP-R34	0.975 (0.967 - 0.982)	0.999 (0.999 - 0.999)	0.926 (0.912 - 0.940)	0.998 (0.997 - 0.998)	0.927 (0.916 - 0.938)	0.925 (0.914 - 0.936)
	Merl-R18	0.973 (0.964 - 0.980)	0.999 (0.999 - 0.999)	0.931 (0.917 - 0.945)	0.998 (0.997 - 0.998)	0.927 (0.917 - 0.938)	0.925 (0.915 - 0.936)
1° AVB	Random Init-R18	0.764 (0.731 - 0.792)	0.994 (0.992 - 0.995)	0.683 (0.655 - 0.711)	0.995 (0.994 - 0.995)	0.715 (0.691 - 0.737)	0.709 (0.685 - 0.732)
	ECGCLIP-R18	0.774 (0.744 - 0.801)	0.994 (0.992 - 0.995)	0.717 (0.688 - 0.744)	0.995 (0.994 - 0.996)	0.738 (0.716 - 0.760)	0.733 (0.710 - 0.755)
	ECGCLIP-R34	0.767 (0.735 - 0.796)	0.994 (0.992 - 0.995)	0.656 (0.627 - 0.684)	0.996 (0.995 - 0.996)	0.710 (0.686 - 0.733)	0.707 (0.682 - 0.730)
	Merl-R18	0.764 (0.732 - 0.793)	0.994 (0.992 - 0.995)	0.730 (0.703 - 0.755)	0.994 (0.993 - 0.995)	0.731 (0.709 - 0.752)	0.725 (0.703 - 0.747)
QT Prolong	Random Init-R18	0.330 (0.299 - 0.363)	0.956 (0.951 - 0.961)	0.132 (0.108 - 0.155)	0.997 (0.997 - 0.998)	0.207 (0.172 - 0.238)	0.246 (0.209 - 0.281)
	ECGCLIP-R18	0.376 (0.341 - 0.411)	0.962 (0.958 - 0.966)	0.169 (0.143 - 0.194)	0.998 (0.997 - 0.998)	0.260 (0.224 - 0.292)	0.301 (0.263 - 0.336)
	ECGCLIP-R34	0.427 (0.391 - 0.465)	0.970 (0.966 - 0.973)	0.255 (0.225 - 0.286)	0.997 (0.996 - 0.997)	0.356 (0.320 - 0.391)	0.380 (0.346 - 0.414)
	Merl-R18	0.338 (0.306 - 0.371)	0.958 (0.953 - 0.963)	0.138 (0.114 - 0.162)	0.997 (0.997 - 0.998)	0.214 (0.182 - 0.245)	0.248 (0.211 - 0.282)
ER	Random Init-R18	0.511 (0.456 - 0.568)	0.984 (0.981 - 0.987)	0.454 (0.405 - 0.508)	0.997 (0.996 - 0.997)	0.496 (0.450 - 0.543)	0.494 (0.448 - 0.541)
	ECGCLIP-R18	0.568 (0.517 - 0.621)	0.989 (0.986 - 0.991)	0.339 (0.292 - 0.387)	0.999 (0.999 - 0.999)	0.463 (0.411 - 0.513)	0.494 (0.446 - 0.543)
	ECGCLIP-R34	0.563 (0.514 - 0.618)	0.991 (0.989 - 0.992)	0.497 (0.448 - 0.551)	0.997 (0.996 - 0.997)	0.533 (0.488 - 0.578)	0.531 (0.487 - 0.577)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAFB	Merl-R18	0.473 (0.414 - 0.531)	0.979 (0.974 - 0.984)	0.484 (0.435 - 0.537)	0.995 (0.995 - 0.996)	0.472 (0.429 - 0.516)	0.468 (0.424 - 0.512)
	Random Init-R18	0.406 (0.307 - 0.518)	0.995 (0.994 - 0.996)	0.510 (0.411 - 0.614)	0.998 (0.998 - 0.999)	0.431 (0.350 - 0.511)	0.435 (0.356 - 0.516)
	ECGCLIP-R18	0.432 (0.329 - 0.543)	0.997 (0.995 - 0.998)	0.500 (0.410 - 0.604)	0.998 (0.998 - 0.999)	0.434 (0.356 - 0.513)	0.436 (0.359 - 0.514)
	ECGCLIP-R34	0.479 (0.380 - 0.579)	0.997 (0.996 - 0.998)	0.510 (0.411 - 0.614)	0.998 (0.998 - 0.999)	0.439 (0.360 - 0.516)	0.442 (0.363 - 0.519)
LBBB	Merl-R18	0.399 (0.300 - 0.511)	0.996 (0.994 - 0.997)	0.500 (0.407 - 0.603)	0.998 (0.997 - 0.998)	0.405 (0.328 - 0.482)	0.411 (0.336 - 0.490)
	Random Init-R18	0.984 (0.970 - 0.994)	1.000 (1.000 - 1.000)	0.971 (0.940 - 0.994)	1.000 (1.000 - 1.000)	0.951 (0.921 - 0.975)	0.951 (0.922 - 0.975)
	ECGCLIP-R18	0.988 (0.975 - 0.996)	1.000 (1.000 - 1.000)	0.957 (0.921 - 0.987)	1.000 (1.000 - 1.000)	0.950 (0.920 - 0.973)	0.950 (0.920 - 0.973)
	ECGCLIP-R34	0.991 (0.981 - 0.998)	1.000 (1.000 - 1.000)	0.964 (0.932 - 0.992)	1.000 (1.000 - 1.000)	0.957 (0.928 - 0.978)	0.957 (0.928 - 0.979)
IVB	Merl-R18	0.985 (0.971 - 0.995)	1.000 (1.000 - 1.000)	0.971 (0.940 - 0.994)	1.000 (1.000 - 1.000)	0.957 (0.930 - 0.978)	0.957 (0.931 - 0.978)
	Random Init-R18	0.700 (0.585 - 0.801)	0.999 (0.998 - 0.999)	0.667 (0.542 - 0.780)	1.000 (0.999 - 1.000)	0.655 (0.554 - 0.745)	0.655 (0.554 - 0.746)
	ECGCLIP-R18	0.759 (0.630 - 0.867)	0.999 (0.998 - 1.000)	0.789 (0.673 - 0.891)	1.000 (0.999 - 1.000)	0.738 (0.645 - 0.814)	0.739 (0.648 - 0.817)
	ECGCLIP-R34	0.801 (0.687 - 0.889)	1.000 (0.999 - 1.000)	0.737 (0.622 - 0.845)	1.000 (1.000 - 1.000)	0.750 (0.660 - 0.827)	0.750 (0.660 - 0.827)
Short PR	Merl-R18	0.685 (0.561 - 0.792)	0.998 (0.997 - 0.999)	0.632 (0.500 - 0.750)	1.000 (0.999 - 1.000)	0.643 (0.530 - 0.736)	0.643 (0.532 - 0.737)
	Random Init-R18	0.133 (0.091 - 0.203)	0.989 (0.986 - 0.993)	0.287 (0.188 - 0.388)	0.997 (0.996 - 0.997)	0.191 (0.124 - 0.259)	0.201 (0.131 - 0.274)
	ECGCLIP-R18	0.169 (0.111 - 0.250)	0.988 (0.982 - 0.993)	0.412 (0.301 - 0.519)	0.997 (0.996 - 0.997)	0.254 (0.187 - 0.326)	0.273 (0.201 - 0.345)
	ECGCLIP-R34	0.199 (0.129 - 0.280)	0.988 (0.978 - 0.994)	0.375 (0.278 - 0.485)	0.998 (0.997 - 0.998)	0.273 (0.196 - 0.349)	0.282 (0.203 - 0.359)
	Merl-R18	0.143 (0.096 - 0.218)	0.988 (0.982 - 0.992)	0.438 (0.329 - 0.548)	0.995 (0.995 - 0.996)	0.219 (0.159 - 0.282)	0.250 (0.188 - 0.316)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
VPE	Random Init-R18	0.826 (0.742 – 0.902)	0.997 (0.994 – 0.999)	0.671 (0.571 – 0.777)	1.000 (1.000 – 1.000)	0.780 (0.701 – 0.855)	0.790 (0.719 – 0.860)
	ECGCLIP-R18	0.829 (0.745 – 0.904)	0.997 (0.994 – 0.999)	0.707 (0.611 – 0.811)	1.000 (1.000 – 1.000)	0.806 (0.732 – 0.877)	0.813 (0.745 – 0.880)
	ECGCLIP-R34	0.826 (0.738 – 0.903)	0.996 (0.993 – 0.999)	0.683 (0.580 – 0.785)	1.000 (1.000 – 1.000)	0.794 (0.718 – 0.864)	0.805 (0.736 – 0.871)
	Merl-R18	0.760 (0.673 – 0.842)	0.993 (0.985 – 0.998)	0.646 (0.544 – 0.753)	1.000 (1.000 – 1.000)	0.741 (0.658 – 0.818)	0.749 (0.672 – 0.823)
3° AVB	Random Init-R18	0.032 (0.013 – 0.084)	0.999 (0.999 – 1.000)	0.000 (0.000 – 0.000)	1.000 (1.000 – 1.000)	0.000 (0.000 – 0.000)	0.000 (0.000 – 0.000)
	ECGCLIP-R18	0.184 (0.033 – 1.000)	1.000 (1.000 – 1.000)	1.000 (1.000 – 1.000)	1.000 (1.000 – 1.000)	0.250 (0.105 – 0.524)	0.378 (0.236 – 0.596)
	ECGCLIP-R34	0.663 (0.100 – 1.000)	1.000 (1.000 – 1.000)	1.000 (1.000 – 1.000)	1.000 (1.000 – 1.000)	0.400 (0.167 – 0.769)	0.500 (0.301 – 0.791)
	Merl-R18	0.012 (0.005 – 0.032)	0.998 (0.998 – 0.999)	0.000 (0.000 – 0.000)	1.000 (1.000 – 1.000)	0.000 (0.000 – 0.000)	0.000 (0.000 – 0.000)
2° 1 Type AVB	Random Init-R18	0.090 (0.018 – 0.283)	0.995 (0.989 – 0.998)	0.182 (0.000 – 0.429)	1.000 (1.000 – 1.000)	0.211 (0.000 – 0.444)	0.213 (0.000 – 0.447)
	ECGCLIP-R18	0.161 (0.036 – 0.458)	0.997 (0.994 – 0.999)	0.364 (0.100 – 0.667)	1.000 (0.999 – 1.000)	0.258 (0.067 – 0.462)	0.269 (0.072 – 0.467)
	ECGCLIP-R34	0.166 (0.044 – 0.373)	0.998 (0.995 – 1.000)	0.273 (0.000 – 0.556)	1.000 (1.000 – 1.000)	0.240 (0.000 – 0.445)	0.242 (0.000 – 0.456)
	Merl-R18	0.158 (0.018 – 0.405)	0.995 (0.990 – 0.998)	0.182 (0.000 – 0.429)	1.000 (1.000 – 1.000)	0.190 (0.000 – 0.429)	0.191 (0.000 – 0.433)
2° 2 Type AVB	Random Init-R18	0.061 (0.010 – 0.197)	0.999 (0.998 – 1.000)	0.400 (0.000 – 1.000)	0.999 (0.999 – 0.999)	0.085 (0.000 – 0.204)	0.138 (0.000 – 0.286)
	ECGCLIP-R18	0.182 (0.018 – 0.679)	0.999 (0.999 – 1.000)	0.200 (0.000 – 0.667)	1.000 (1.000 – 1.000)	0.222 (0.000 – 0.571)	0.224 (0.000 – 0.577)
	ECGCLIP-R34	0.343 (0.070 – 0.944)	1.000 (0.999 – 1.000)	0.600 (0.000 – 1.000)	1.000 (1.000 – 1.000)	0.500 (0.000 – 0.800)	0.507 (0.000 – 0.816)
	Merl-R18	0.096 (0.013 – 0.546)	0.999 (0.998 – 1.000)	0.200 (0.000 – 0.667)	1.000 (1.000 – 1.000)	0.222 (0.000 – 0.600)	0.224 (0.000 – 0.612)
LVH	Random Init-R18	0.847 (0.833 – 0.861)	0.990 (0.989 – 0.991)	0.509 (0.487 – 0.531)	0.998 (0.998 – 0.999)	0.656 (0.635 – 0.677)	0.676 (0.658 – 0.694)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAH	ECGCLIP-R18	0.870 (0.858 - 0.882)	0.992 (0.991 - 0.993)	0.550 (0.528 - 0.572)	0.998 (0.998 - 0.999)	0.690 (0.671 - 0.709)	0.706 (0.688 - 0.722)
	ECGCLIP-R34	0.863 (0.850 - 0.877)	0.992 (0.992 - 0.993)	0.549 (0.527 - 0.571)	0.998 (0.997 - 0.998)	0.684 (0.667 - 0.703)	0.698 (0.681 - 0.715)
	Merl-R18	0.834 (0.820 - 0.848)	0.989 (0.988 - 0.990)	0.524 (0.503 - 0.546)	0.998 (0.997 - 0.998)	0.662 (0.643 - 0.682)	0.677 (0.660 - 0.694)
	Random Init-R18	0.668 (0.642 - 0.695)	0.977 (0.974 - 0.981)	0.297 (0.273 - 0.324)	0.998 (0.998 - 0.999)	0.438 (0.409 - 0.468)	0.490 (0.465 - 0.517)
	ECGCLIP-R18	0.673 (0.647 - 0.701)	0.980 (0.978 - 0.983)	0.291 (0.266 - 0.317)	0.998 (0.998 - 0.999)	0.431 (0.401 - 0.462)	0.485 (0.457 - 0.512)
	ECGCLIP-R34	0.681 (0.656 - 0.709)	0.982 (0.980 - 0.984)	0.353 (0.328 - 0.381)	0.997 (0.997 - 0.998)	0.490 (0.463 - 0.519)	0.523 (0.498 - 0.551)
	Merl-R18	0.676 (0.648 - 0.702)	0.978 (0.974 - 0.981)	0.346 (0.320 - 0.373)	0.998 (0.997 - 0.998)	0.484 (0.456 - 0.514)	0.520 (0.495 - 0.547)
	Random Init-R18	0.423 (0.357 - 0.498)	0.991 (0.987 - 0.994)	0.411 (0.344 - 0.476)	0.998 (0.998 - 0.998)	0.443 (0.382 - 0.504)	0.442 (0.381 - 0.504)
	ECGCLIP-R18	0.527 (0.456 - 0.603)	0.994 (0.991 - 0.997)	0.484 (0.416 - 0.560)	0.998 (0.998 - 0.999)	0.520 (0.458 - 0.580)	0.519 (0.456 - 0.579)
	ECGCLIP-R34	0.610 (0.540 - 0.680)	0.996 (0.994 - 0.998)	0.526 (0.456 - 0.596)	0.998 (0.998 - 0.999)	0.547 (0.484 - 0.607)	0.546 (0.483 - 0.606)
	Merl-R18	0.459 (0.387 - 0.534)	0.993 (0.991 - 0.995)	0.339 (0.275 - 0.406)	0.999 (0.999 - 0.999)	0.428 (0.357 - 0.497)	0.441 (0.370 - 0.509)
	Random Init-R18	0.226 (0.092 - 0.416)	0.999 (0.998 - 0.999)	0.421 (0.190 - 0.647)	0.999 (0.999 - 0.999)	0.211 (0.083 - 0.333)	0.243 (0.105 - 0.373)
VA	ECGCLIP-R18	0.258 (0.100 - 0.467)	0.999 (0.998 - 0.999)	0.421 (0.200 - 0.643)	0.999 (0.999 - 0.999)	0.219 (0.095 - 0.353)	0.249 (0.112 - 0.383)
	ECGCLIP-R34	0.354 (0.155 - 0.573)	0.999 (0.998 - 0.999)	0.474 (0.250 - 0.700)	0.999 (0.999 - 0.999)	0.265 (0.127 - 0.410)	0.295 (0.146 - 0.440)
	Merl-R18	0.124 (0.064 - 0.226)	0.998 (0.997 - 0.999)	0.421 (0.200 - 0.643)	0.999 (0.998 - 0.999)	0.190 (0.078 - 0.304)	0.227 (0.101 - 0.356)
RAH	Random Init-R18	0.130 (0.051 - 0.237)	0.976 (0.968 - 0.985)	0.163 (0.068 - 0.264)	1.000 (0.999 - 1.000)	0.205 (0.087 - 0.320)	0.212 (0.089 - 0.335)
	ECGCLIP-R18	0.155 (0.070 - 0.268)	0.985 (0.980 - 0.991)	0.163 (0.068 - 0.273)	1.000 (0.999 - 1.000)	0.211 (0.088 - 0.337)	0.219 (0.094 - 0.351)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
DC	ECGCLIP-R34	0.234 (0.127 - 0.367)	0.991 (0.987 - 0.994)	0.163 (0.068 - 0.269)	1.000 (1.000 - 1.000)	0.262 (0.118 - 0.406)	0.330 (0.173 - 0.473)
	Merl-R18	0.147 (0.063 - 0.253)	0.980 (0.974 - 0.986)	0.082 (0.019 - 0.164)	1.000 (1.000 - 1.000)	0.140 (0.034 - 0.262)	0.202 (0.054 - 0.342)
	Random Init-R18	0.000 (0.000 - 0.001)	0.627 (0.414 - 0.811)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.000 (0.000 - 0.001)	0.710 (0.561 - 0.855)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.000 (0.000 - 0.001)	0.783 (0.685 - 0.866)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	Merl-R18	0.000 (0.000 - 0.001)	0.759 (0.580 - 0.892)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.001)	0.000 (0.000 - 0.000)
VVI	Random Init-R18	0.515 (0.358 - 0.665)	0.998 (0.996 - 0.999)	0.500 (0.348 - 0.650)	1.000 (0.999 - 1.000)	0.512 (0.364 - 0.631)	0.511 (0.363 - 0.632)
	ECGCLIP-R18	0.567 (0.421 - 0.703)	0.998 (0.996 - 0.999)	0.409 (0.275 - 0.563)	1.000 (1.000 - 1.000)	0.480 (0.338 - 0.620)	0.487 (0.343 - 0.628)
	ECGCLIP-R34	0.653 (0.517 - 0.771)	0.998 (0.997 - 0.999)	0.568 (0.421 - 0.718)	1.000 (1.000 - 1.000)	0.610 (0.475 - 0.732)	0.611 (0.477 - 0.731)
	Merl-R18	0.491 (0.340 - 0.635)	0.997 (0.996 - 0.999)	0.295 (0.171 - 0.432)	1.000 (1.000 - 1.000)	0.394 (0.241 - 0.533)	0.417 (0.273 - 0.555)
VAT	Random Init-R18	0.675 (0.477 - 0.833)	0.956 (0.902 - 0.996)	0.625 (0.454 - 0.786)	1.000 (1.000 - 1.000)	0.702 (0.550 - 0.821)	0.707 (0.567 - 0.823)
	ECGCLIP-R18	0.761 (0.590 - 0.892)	0.977 (0.946 - 0.997)	0.719 (0.556 - 0.870)	1.000 (1.000 - 1.000)	0.793 (0.666 - 0.897)	0.797 (0.669 - 0.899)
	ECGCLIP-R34	0.794 (0.644 - 0.911)	0.994 (0.987 - 0.999)	0.688 (0.517 - 0.840)	1.000 (1.000 - 1.000)	0.800 (0.666 - 0.897)	0.811 (0.691 - 0.901)
DDD	Merl-R18	0.704 (0.522 - 0.848)	0.993 (0.986 - 0.998)	0.625 (0.455 - 0.794)	1.000 (1.000 - 1.000)	0.702 (0.549 - 0.824)	0.707 (0.556 - 0.825)
	Random Init-R18	0.122 (0.030 - 0.319)	1.000 (1.000 - 1.000)	0.667 (0.000 - 1.000)	1.000 (1.000 - 1.000)	0.267 (0.000 - 0.571)	0.333 (0.000 - 0.608)
	ECGCLIP-R18	0.261 (0.028 - 1.000)	1.000 (1.000 - 1.000)	0.667 (0.000 - 1.000)	1.000 (1.000 - 1.000)	0.222 (0.000 - 0.476)	0.298 (0.000 - 0.539)
	ECGCLIP-R34	0.288 (0.056 - 1.000)	1.000 (1.000 - 1.000)	0.667 (0.000 - 1.000)	1.000 (1.000 - 1.000)	0.286 (0.000 - 0.571)	0.348 (0.000 - 0.632)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AAI	Merl-R18	0.313 (0.024 - 1.000)	1.000 (0.999 - 1.000)	0.667 (0.000 - 1.000)	1.000 (1.000 - 1.000)	0.444 (0.000 - 0.800)	0.471 (0.000 - 0.816)
	Random Init-R18	0.189 (0.070 - 0.342)	0.980 (0.955 - 0.998)	0.200 (0.050 - 0.360)	1.000 (1.000 - 1.000)	0.238 (0.061 - 0.400)	0.242 (0.067 - 0.411)
	ECGCLIP-R18	0.423 (0.208 - 0.620)	0.989 (0.976 - 0.999)	0.280 (0.100 - 0.458)	1.000 (1.000 - 1.000)	0.424 (0.176 - 0.621)	0.495 (0.288 - 0.667)
	ECGCLIP-R34	0.781 (0.614 - 0.922)	0.993 (0.981 - 1.000)	0.640 (0.440 - 0.815)	1.000 (1.000 - 1.000)	0.762 (0.588 - 0.889)	0.776 (0.628 - 0.891)
PSM	Merl-R18	0.271 (0.093 - 0.455)	0.985 (0.968 - 0.998)	0.160 (0.031 - 0.313)	1.000 (1.000 - 1.000)	0.267 (0.060 - 0.457)	0.358 (0.104 - 0.534)
	Random Init-R18	0.008 (0.000 - 0.026)	0.984 (0.956 - 0.998)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.005 (0.000 - 0.017)	0.957 (0.876 - 0.998)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.003 (0.000 - 0.011)	0.980 (0.950 - 0.997)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
HK	Merl-R18	0.005 (0.000 - 0.025)	0.970 (0.918 - 0.999)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Random Init-R18	0.202 (0.000 - 0.667)	0.975 (0.954 - 0.993)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.202 (0.000 - 0.667)	0.975 (0.955 - 0.992)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.026 (0.002 - 0.122)	0.993 (0.985 - 0.999)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Merl-R18	0.022 (0.001 - 0.168)	0.976 (0.958 - 0.992)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S20: Detailed task-specific performance comparison on the external MIMIC-IV-ECG cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
OMI	Random Init-R18	0.060 (0.059 – 0.062)	0.795 (0.792 – 0.799)	0.074 (0.070 – 0.079)	0.982 (0.982 – 0.982)	0.074 (0.070 – 0.079)	0.056 (0.052 – 0.061)
	ECGCLIP-R18	0.070 (0.068 – 0.071)	0.817 (0.814 – 0.820)	0.084 (0.079 – 0.088)	0.984 (0.983 – 0.984)	0.087 (0.082 – 0.092)	0.070 (0.065 – 0.075)
	ECGCLIP-R34	0.069 (0.067 – 0.071)	0.823 (0.820 – 0.826)	0.074 (0.070 – 0.079)	0.984 (0.984 – 0.984)	0.078 (0.074 – 0.083)	0.062 (0.057 – 0.066)
	Merl-R18	0.152 (0.147 – 0.158)	0.846 (0.843 – 0.849)	0.070 (0.066 – 0.074)	0.998 (0.998 – 0.998)	0.119 (0.112 – 0.125)	0.159 (0.150 – 0.167)
STEMI	Random Init-R18	0.578 (0.575 – 0.580)	0.807 (0.806 – 0.808)	0.019 (0.019 – 0.020)	0.999 (0.999 – 0.999)	0.038 (0.037 – 0.039)	0.112 (0.110 – 0.115)
	ECGCLIP-R18	0.596 (0.594 – 0.599)	0.821 (0.820 – 0.822)	0.016 (0.016 – 0.017)	1.000 (1.000 – 1.000)	0.032 (0.031 – 0.033)	0.106 (0.104 – 0.108)
	ECGCLIP-R34	0.585 (0.582 – 0.588)	0.820 (0.819 – 0.821)	0.026 (0.026 – 0.027)	0.999 (0.999 – 0.999)	0.051 (0.050 – 0.053)	0.133 (0.130 – 0.135)
	Merl-R18	0.575 (0.573 – 0.578)	0.823 (0.822 – 0.825)	0.012 (0.011 – 0.012)	0.999 (0.999 – 0.999)	0.023 (0.022 – 0.024)	0.075 (0.072 – 0.078)
NSR	Random Init-R18	0.367 (0.365 – 0.370)	0.796 (0.794 – 0.797)	0.978 (0.977 – 0.979)	0.420 (0.419 – 0.421)	0.378 (0.377 – 0.380)	0.299 (0.298 – 0.300)
	ECGCLIP-R18	0.345 (0.343 – 0.348)	0.789 (0.788 – 0.790)	0.963 (0.962 – 0.964)	0.449 (0.448 – 0.451)	0.386 (0.384 – 0.387)	0.306 (0.304 – 0.307)
	ECGCLIP-R34	0.344 (0.342 – 0.346)	0.789 (0.787 – 0.790)	0.980 (0.979 – 0.980)	0.429 (0.428 – 0.430)	0.382 (0.381 – 0.384)	0.306 (0.305 – 0.307)
	Merl-R18	0.325 (0.323 – 0.327)	0.787 (0.786 – 0.788)	0.945 (0.943 – 0.946)	0.471 (0.470 – 0.473)	0.389 (0.388 – 0.391)	0.305 (0.304 – 0.307)
SBrad	Random Init-R18	0.868 (0.865 – 0.870)	0.980 (0.980 – 0.981)	0.835 (0.832 – 0.837)	0.970 (0.970 – 0.970)	0.819 (0.817 – 0.821)	0.792 (0.790 – 0.794)
	ECGCLIP-R18	0.868 (0.866 – 0.870)	0.981 (0.981 – 0.981)	0.810 (0.808 – 0.813)	0.974 (0.973 – 0.974)	0.815 (0.813 – 0.817)	0.788 (0.786 – 0.790)
	ECGCLIP-R34	0.891 (0.889 – 0.893)	0.987 (0.986 – 0.987)	0.874 (0.872 – 0.876)	0.978 (0.977 – 0.978)	0.862 (0.861 – 0.864)	0.842 (0.840 – 0.844)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STach	Merl-R18	0.841 (0.838 - 0.843)	0.976 (0.976 - 0.977)	0.638 (0.635 - 0.641)	0.984 (0.983 - 0.984)	0.729 (0.727 - 0.732)	0.705 (0.703 - 0.708)
	Random Init-R18	0.917 (0.915 - 0.919)	0.988 (0.988 - 0.989)	0.795 (0.792 - 0.798)	0.993 (0.993 - 0.994)	0.853 (0.851 - 0.855)	0.843 (0.841 - 0.845)
	ECGCLIP-R18	0.920 (0.918 - 0.922)	0.987 (0.986 - 0.987)	0.801 (0.798 - 0.803)	0.994 (0.994 - 0.994)	0.859 (0.857 - 0.861)	0.849 (0.847 - 0.851)
	ECGCLIP-R34	0.937 (0.936 - 0.939)	0.991 (0.990 - 0.991)	0.831 (0.828 - 0.834)	0.995 (0.994 - 0.995)	0.881 (0.879 - 0.882)	0.872 (0.870 - 0.874)
VPB	Merl-R18	0.911 (0.910 - 0.913)	0.986 (0.986 - 0.987)	0.761 (0.758 - 0.764)	0.994 (0.993 - 0.994)	0.833 (0.830 - 0.835)	0.822 (0.820 - 0.825)
	Random Init-R18	0.554 (0.549 - 0.558)	0.912 (0.910 - 0.913)	0.550 (0.546 - 0.554)	0.971 (0.971 - 0.971)	0.562 (0.559 - 0.566)	0.532 (0.528 - 0.536)
	ECGCLIP-R18	0.608 (0.604 - 0.612)	0.925 (0.924 - 0.926)	0.577 (0.573 - 0.582)	0.977 (0.976 - 0.977)	0.607 (0.604 - 0.611)	0.581 (0.578 - 0.585)
	ECGCLIP-R34	0.713 (0.709 - 0.717)	0.947 (0.946 - 0.948)	0.690 (0.686 - 0.694)	0.981 (0.981 - 0.982)	0.708 (0.705 - 0.711)	0.688 (0.685 - 0.691)
APB	Merl-R18	0.518 (0.514 - 0.523)	0.897 (0.895 - 0.898)	0.540 (0.536 - 0.544)	0.965 (0.965 - 0.965)	0.533 (0.529 - 0.536)	0.499 (0.495 - 0.502)
	Random Init-R18	0.135 (0.132 - 0.137)	0.761 (0.758 - 0.763)	0.614 (0.609 - 0.619)	0.768 (0.767 - 0.769)	0.191 (0.188 - 0.193)	0.184 (0.182 - 0.187)
	ECGCLIP-R18	0.181 (0.178 - 0.185)	0.787 (0.784 - 0.789)	0.577 (0.572 - 0.582)	0.829 (0.828 - 0.830)	0.225 (0.222 - 0.228)	0.217 (0.214 - 0.220)
	ECGCLIP-R34	0.330 (0.325 - 0.335)	0.845 (0.842 - 0.847)	0.565 (0.560 - 0.571)	0.922 (0.921 - 0.923)	0.355 (0.351 - 0.358)	0.339 (0.335 - 0.343)
SArr	Merl-R18	0.135 (0.132 - 0.137)	0.760 (0.758 - 0.763)	0.682 (0.677 - 0.687)	0.713 (0.712 - 0.714)	0.178 (0.176 - 0.180)	0.179 (0.177 - 0.181)
	Random Init-R18	0.053 (0.052 - 0.054)	0.625 (0.622 - 0.628)	0.021 (0.019 - 0.023)	0.993 (0.993 - 0.993)	0.034 (0.031 - 0.037)	0.029 (0.025 - 0.033)
	ECGCLIP-R18	0.054 (0.053 - 0.055)	0.631 (0.628 - 0.634)	0.028 (0.026 - 0.030)	0.991 (0.990 - 0.991)	0.043 (0.040 - 0.046)	0.034 (0.030 - 0.038)
	ECGCLIP-R34	0.055 (0.054 - 0.056)	0.638 (0.635 - 0.642)	0.069 (0.066 - 0.072)	0.973 (0.972 - 0.973)	0.074 (0.071 - 0.078)	0.045 (0.042 - 0.048)
	Merl-R18	0.055 (0.054 - 0.056)	0.633 (0.630 - 0.637)	0.097 (0.093 - 0.101)	0.960 (0.960 - 0.961)	0.087 (0.084 - 0.090)	0.052 (0.049 - 0.055)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AF	Random Init-R18	0.715 (0.711 - 0.719)	0.955 (0.954 - 0.956)	0.496 (0.493 - 0.500)	0.986 (0.986 - 0.987)	0.615 (0.612 - 0.618)	0.603 (0.599 - 0.606)
	ECGCLIP-R18	0.724 (0.720 - 0.728)	0.953 (0.952 - 0.954)	0.463 (0.459 - 0.467)	0.989 (0.988 - 0.989)	0.593 (0.589 - 0.596)	0.588 (0.585 - 0.591)
	ECGCLIP-R34	0.725 (0.721 - 0.728)	0.948 (0.947 - 0.949)	0.489 (0.485 - 0.493)	0.988 (0.987 - 0.988)	0.613 (0.609 - 0.616)	0.603 (0.600 - 0.607)
	Merl-R18	0.713 (0.709 - 0.717)	0.953 (0.952 - 0.954)	0.510 (0.506 - 0.513)	0.985 (0.985 - 0.985)	0.622 (0.619 - 0.625)	0.606 (0.603 - 0.610)
JPB	Random Init-R18	0.006 (0.006 - 0.007)	0.691 (0.680 - 0.702)	0.216 (0.198 - 0.234)	0.930 (0.930 - 0.931)	0.017 (0.016 - 0.019)	0.031 (0.027 - 0.034)
	ECGCLIP-R18	0.008 (0.007 - 0.008)	0.720 (0.709 - 0.730)	0.131 (0.116 - 0.145)	0.971 (0.971 - 0.972)	0.024 (0.021 - 0.027)	0.033 (0.028 - 0.037)
	ECGCLIP-R34	0.014 (0.013 - 0.016)	0.770 (0.760 - 0.781)	0.156 (0.141 - 0.172)	0.985 (0.985 - 0.985)	0.050 (0.045 - 0.055)	0.062 (0.055 - 0.069)
	Merl-R18	0.007 (0.006 - 0.007)	0.693 (0.682 - 0.703)	0.038 (0.030 - 0.047)	0.992 (0.992 - 0.992)	0.020 (0.015 - 0.024)	0.018 (0.013 - 0.023)
AJR	Random Init-R18	0.012 (0.011 - 0.013)	0.743 (0.735 - 0.751)	0.109 (0.098 - 0.120)	0.982 (0.982 - 0.983)	0.040 (0.036 - 0.044)	0.043 (0.038 - 0.049)
	ECGCLIP-R18	0.013 (0.012 - 0.014)	0.731 (0.723 - 0.739)	0.118 (0.107 - 0.129)	0.982 (0.982 - 0.982)	0.042 (0.038 - 0.046)	0.047 (0.042 - 0.052)
	ECGCLIP-R34	0.013 (0.012 - 0.014)	0.756 (0.748 - 0.764)	0.077 (0.067 - 0.086)	0.989 (0.989 - 0.989)	0.040 (0.035 - 0.045)	0.039 (0.033 - 0.045)
	Merl-R18	0.011 (0.010 - 0.012)	0.757 (0.749 - 0.764)	0.031 (0.025 - 0.038)	0.992 (0.992 - 0.992)	0.020 (0.017 - 0.025)	0.016 (0.012 - 0.021)
ATach	Random Init-R18	0.035 (0.031 - 0.040)	0.857 (0.848 - 0.866)	0.416 (0.394 - 0.439)	0.967 (0.966 - 0.967)	0.056 (0.052 - 0.060)	0.104 (0.098 - 0.110)
	ECGCLIP-R18	0.031 (0.028 - 0.035)	0.848 (0.839 - 0.857)	0.330 (0.308 - 0.350)	0.981 (0.980 - 0.981)	0.072 (0.067 - 0.077)	0.110 (0.102 - 0.117)
	ECGCLIP-R34	0.024 (0.021 - 0.026)	0.891 (0.885 - 0.897)	0.235 (0.216 - 0.252)	0.984 (0.983 - 0.984)	0.060 (0.055 - 0.065)	0.084 (0.077 - 0.091)
	Merl-R18	0.025 (0.022 - 0.028)	0.850 (0.842 - 0.858)	0.314 (0.294 - 0.336)	0.976 (0.976 - 0.977)	0.058 (0.053 - 0.062)	0.093 (0.087 - 0.101)
JTach	Random Init-R18	0.004 (0.003 - 0.005)	0.944 (0.931 - 0.955)	0.149 (0.099 - 0.205)	0.995 (0.995 - 0.996)	0.013 (0.009 - 0.019)	0.031 (0.020 - 0.044)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
SVT	ECGCLIP-R18	0.006 (0.005 - 0.008)	0.945 (0.931 - 0.958)	0.238 (0.173 - 0.304)	0.996 (0.996 - 0.996)	0.023 (0.016 - 0.031)	0.053 (0.038 - 0.068)
	ECGCLIP-R34	0.009 (0.006 - 0.012)	0.953 (0.935 - 0.966)	0.327 (0.256 - 0.397)	0.996 (0.996 - 0.996)	0.031 (0.023 - 0.039)	0.072 (0.056 - 0.089)
	Merl-R18	0.012 (0.008 - 0.018)	0.935 (0.918 - 0.951)	0.196 (0.138 - 0.257)	0.998 (0.998 - 0.998)	0.038 (0.026 - 0.052)	0.064 (0.044 - 0.086)
	Random Init-R18	0.402 (0.369 - 0.436)	0.987 (0.981 - 0.991)	0.315 (0.286 - 0.343)	1.000 (1.000 - 1.000)	0.404 (0.371 - 0.435)	0.420 (0.389 - 0.451)
	ECGCLIP-R18	0.399 (0.364 - 0.430)	0.989 (0.986 - 0.992)	0.299 (0.271 - 0.327)	1.000 (1.000 - 1.000)	0.402 (0.370 - 0.433)	0.428 (0.396 - 0.459)
	ECGCLIP-R34	0.380 (0.347 - 0.415)	0.991 (0.987 - 0.993)	0.290 (0.261 - 0.319)	1.000 (1.000 - 1.000)	0.383 (0.352 - 0.415)	0.405 (0.373 - 0.437)
	Merl-R18	0.342 (0.308 - 0.373)	0.980 (0.975 - 0.985)	0.180 (0.156 - 0.202)	1.000 (1.000 - 1.000)	0.286 (0.254 - 0.317)	0.354 (0.321 - 0.386)
	Random Init-R18	0.005 (0.004 - 0.006)	0.872 (0.861 - 0.883)	0.021 (0.009 - 0.033)	0.998 (0.998 - 0.998)	0.011 (0.005 - 0.018)	0.011 (0.004 - 0.019)
VEB	ECGCLIP-R18	0.005 (0.004 - 0.006)	0.874 (0.862 - 0.885)	0.007 (0.002 - 0.014)	0.999 (0.999 - 0.999)	0.006 (0.001 - 0.012)	0.005 (0.001 - 0.012)
	ECGCLIP-R34	0.004 (0.004 - 0.005)	0.871 (0.860 - 0.880)	0.005 (0.000 - 0.012)	0.999 (0.999 - 0.999)	0.004 (0.000 - 0.009)	0.003 (-0.001 - 0.008)
	Merl-R18	0.005 (0.004 - 0.006)	0.873 (0.862 - 0.883)	0.067 (0.047 - 0.090)	0.994 (0.994 - 0.994)	0.016 (0.011 - 0.021)	0.022 (0.015 - 0.030)
	Random Init-R18	0.046 (0.038 - 0.056)	0.952 (0.947 - 0.958)	0.201 (0.174 - 0.230)	0.998 (0.998 - 0.998)	0.126 (0.107 - 0.144)	0.134 (0.115 - 0.154)
VTach	ECGCLIP-R18	0.049 (0.040 - 0.060)	0.956 (0.949 - 0.961)	0.190 (0.161 - 0.219)	0.998 (0.998 - 0.998)	0.126 (0.108 - 0.146)	0.133 (0.113 - 0.154)
	ECGCLIP-R34	0.053 (0.045 - 0.064)	0.963 (0.957 - 0.968)	0.244 (0.215 - 0.277)	0.998 (0.998 - 0.998)	0.134 (0.117 - 0.152)	0.149 (0.131 - 0.169)
	Merl-R18	0.049 (0.041 - 0.061)	0.948 (0.940 - 0.955)	0.101 (0.078 - 0.124)	0.999 (0.999 - 0.999)	0.111 (0.087 - 0.136)	0.111 (0.086 - 0.136)
RBBB	Random Init-R18	0.900 (0.898 - 0.902)	0.985 (0.985 - 0.985)	0.875 (0.873 - 0.878)	0.973 (0.973 - 0.974)	0.807 (0.804 - 0.809)	0.790 (0.788 - 0.793)
	ECGCLIP-R18	0.898 (0.896 - 0.900)	0.984 (0.984 - 0.985)	0.852 (0.850 - 0.855)	0.978 (0.977 - 0.978)	0.813 (0.811 - 0.815)	0.796 (0.794 - 0.798)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
1° AVB	ECGCLIP-R34	0.902 (0.900 - 0.904)	0.984 (0.984 - 0.984)	0.866 (0.864 - 0.869)	0.973 (0.973 - 0.974)	0.802 (0.800 - 0.804)	0.785 (0.782 - 0.787)
	Merl-R18	0.891 (0.889 - 0.893)	0.984 (0.984 - 0.985)	0.779 (0.776 - 0.783)	0.989 (0.989 - 0.989)	0.821 (0.818 - 0.823)	0.807 (0.804 - 0.809)
	Random Init-R18	0.648 (0.644 - 0.652)	0.960 (0.959 - 0.961)	0.640 (0.636 - 0.644)	0.971 (0.970 - 0.971)	0.640 (0.637 - 0.643)	0.610 (0.607 - 0.614)
	ECGCLIP-R18	0.643 (0.639 - 0.647)	0.959 (0.958 - 0.960)	0.662 (0.658 - 0.665)	0.969 (0.969 - 0.969)	0.648 (0.645 - 0.651)	0.619 (0.616 - 0.622)
	ECGCLIP-R34	0.650 (0.646 - 0.654)	0.959 (0.958 - 0.960)	0.619 (0.615 - 0.622)	0.972 (0.971 - 0.972)	0.630 (0.627 - 0.633)	0.600 (0.597 - 0.604)
	Merl-R18	0.614 (0.610 - 0.619)	0.958 (0.957 - 0.958)	0.654 (0.650 - 0.658)	0.969 (0.968 - 0.969)	0.641 (0.638 - 0.645)	0.612 (0.608 - 0.615)
QT Prolong	Random Init-R18	0.292 (0.288 - 0.296)	0.879 (0.877 - 0.880)	0.194 (0.190 - 0.198)	0.986 (0.986 - 0.986)	0.265 (0.261 - 0.270)	0.261 (0.257 - 0.266)
	ECGCLIP-R18	0.313 (0.309 - 0.318)	0.885 (0.884 - 0.887)	0.211 (0.207 - 0.215)	0.986 (0.986 - 0.986)	0.286 (0.282 - 0.291)	0.283 (0.278 - 0.288)
	ECGCLIP-R34	0.343 (0.338 - 0.348)	0.895 (0.894 - 0.896)	0.286 (0.281 - 0.290)	0.981 (0.981 - 0.981)	0.347 (0.342 - 0.352)	0.328 (0.324 - 0.333)
	Merl-R18	0.243 (0.239 - 0.247)	0.852 (0.851 - 0.854)	0.063 (0.060 - 0.066)	0.996 (0.996 - 0.996)	0.110 (0.106 - 0.115)	0.154 (0.149 - 0.160)
ER	Random Init-R18	0.037 (0.034 - 0.041)	0.722 (0.713 - 0.732)	0.049 (0.043 - 0.056)	0.998 (0.998 - 0.998)	0.070 (0.061 - 0.080)	0.074 (0.064 - 0.085)
	ECGCLIP-R18	0.041 (0.038 - 0.045)	0.751 (0.742 - 0.760)	0.020 (0.016 - 0.025)	0.999 (0.999 - 0.999)	0.034 (0.027 - 0.041)	0.045 (0.035 - 0.054)
	ECGCLIP-R34	0.044 (0.041 - 0.048)	0.795 (0.787 - 0.803)	0.061 (0.054 - 0.068)	0.997 (0.997 - 0.997)	0.077 (0.069 - 0.086)	0.076 (0.067 - 0.085)
LAFB	Merl-R18	0.023 (0.021 - 0.025)	0.707 (0.698 - 0.716)	0.001 (0.000 - 0.002)	1.000 (1.000 - 1.000)	0.002 (0.000 - 0.004)	0.004 (0.000 - 0.008)
	Random Init-R18	0.672 (0.668 - 0.677)	0.961 (0.961 - 0.962)	0.387 (0.383 - 0.391)	0.992 (0.992 - 0.993)	0.515 (0.511 - 0.520)	0.527 (0.522 - 0.531)
	ECGCLIP-R18	0.688 (0.683 - 0.692)	0.962 (0.962 - 0.963)	0.262 (0.258 - 0.266)	0.997 (0.997 - 0.997)	0.402 (0.397 - 0.407)	0.459 (0.455 - 0.463)
	ECGCLIP-R34	0.694 (0.689 - 0.698)	0.963 (0.963 - 0.964)	0.285 (0.281 - 0.289)	0.997 (0.997 - 0.997)	0.427 (0.423 - 0.432)	0.477 (0.473 - 0.481)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LBBB	Merl-R18	0.646 (0.641 - 0.651)	0.950 (0.949 - 0.950)	0.201 (0.197 - 0.205)	0.998 (0.998 - 0.998)	0.325 (0.321 - 0.331)	0.400 (0.395 - 0.404)
	Random Init-R18	0.840 (0.837 - 0.844)	0.980 (0.979 - 0.981)	0.480 (0.474 - 0.485)	0.999 (0.999 - 0.999)	0.637 (0.631 - 0.642)	0.666 (0.662 - 0.671)
	ECGCLIP-R18	0.833 (0.829 - 0.837)	0.979 (0.978 - 0.980)	0.449 (0.443 - 0.455)	0.999 (0.999 - 0.999)	0.610 (0.605 - 0.616)	0.646 (0.641 - 0.651)
	ECGCLIP-R34	0.824 (0.820 - 0.828)	0.977 (0.976 - 0.977)	0.473 (0.467 - 0.478)	0.999 (0.999 - 0.999)	0.629 (0.624 - 0.634)	0.659 (0.654 - 0.663)
IVB	Merl-R18	0.821 (0.817 - 0.825)	0.978 (0.977 - 0.978)	0.249 (0.244 - 0.254)	1.000 (1.000 - 1.000)	0.396 (0.390 - 0.402)	0.484 (0.478 - 0.489)
	Random Init-R18	0.359 (0.354 - 0.365)	0.913 (0.911 - 0.915)	0.261 (0.256 - 0.266)	0.989 (0.988 - 0.989)	0.342 (0.337 - 0.348)	0.341 (0.336 - 0.347)
	ECGCLIP-R18	0.360 (0.355 - 0.366)	0.912 (0.911 - 0.914)	0.243 (0.238 - 0.248)	0.989 (0.989 - 0.990)	0.327 (0.321 - 0.332)	0.330 (0.324 - 0.335)
	ECGCLIP-R34	0.353 (0.348 - 0.359)	0.904 (0.902 - 0.906)	0.225 (0.220 - 0.229)	0.991 (0.991 - 0.991)	0.313 (0.307 - 0.318)	0.322 (0.316 - 0.328)
Short PR	Merl-R18	0.326 (0.321 - 0.331)	0.916 (0.914 - 0.917)	0.120 (0.116 - 0.123)	0.993 (0.993 - 0.993)	0.186 (0.181 - 0.191)	0.209 (0.203 - 0.214)
	Random Init-R18	0.148 (0.142 - 0.155)	0.786 (0.781 - 0.791)	0.098 (0.092 - 0.104)	0.998 (0.998 - 0.998)	0.156 (0.148 - 0.165)	0.189 (0.179 - 0.199)
	ECGCLIP-R18	0.153 (0.146 - 0.160)	0.800 (0.795 - 0.805)	0.071 (0.066 - 0.076)	0.999 (0.999 - 0.999)	0.122 (0.114 - 0.130)	0.170 (0.160 - 0.180)
	ECGCLIP-R34	0.183 (0.175 - 0.190)	0.834 (0.830 - 0.839)	0.095 (0.090 - 0.101)	0.998 (0.998 - 0.999)	0.159 (0.150 - 0.167)	0.207 (0.197 - 0.217)
VPE	Merl-R18	0.100 (0.096 - 0.104)	0.788 (0.783 - 0.792)	0.019 (0.017 - 0.022)	0.999 (0.999 - 0.999)	0.035 (0.030 - 0.039)	0.056 (0.049 - 0.064)
	Random Init-R18	0.133 (0.103 - 0.168)	0.899 (0.882 - 0.914)	0.248 (0.210 - 0.289)	1.000 (1.000 - 1.000)	0.296 (0.256 - 0.336)	0.301 (0.261 - 0.342)
	ECGCLIP-R18	0.145 (0.114 - 0.182)	0.881 (0.864 - 0.896)	0.248 (0.212 - 0.290)	1.000 (1.000 - 1.000)	0.290 (0.252 - 0.331)	0.293 (0.255 - 0.335)
	ECGCLIP-R34	0.141 (0.108 - 0.178)	0.874 (0.856 - 0.891)	0.228 (0.193 - 0.267)	1.000 (1.000 - 1.000)	0.276 (0.237 - 0.317)	0.281 (0.243 - 0.323)
	Merl-R18	0.096 (0.073 - 0.128)	0.870 (0.853 - 0.888)	0.172 (0.141 - 0.208)	1.000 (1.000 - 1.000)	0.239 (0.199 - 0.281)	0.259 (0.218 - 0.302)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
3° AVB	Random Init-R18	0.008 (0.006 – 0.011)	0.940 (0.923 – 0.954)	0.314 (0.250 – 0.381)	0.993 (0.993 – 0.993)	0.023 (0.018 – 0.029)	0.060 (0.047 – 0.074)
	ECGCLIP-R18	0.013 (0.008 – 0.021)	0.939 (0.922 – 0.954)	0.402 (0.338 – 0.474)	0.991 (0.991 – 0.991)	0.023 (0.018 – 0.029)	0.068 (0.056 – 0.081)
	ECGCLIP-R34	0.023 (0.015 – 0.042)	0.948 (0.931 – 0.962)	0.245 (0.186 – 0.310)	0.998 (0.998 – 0.998)	0.063 (0.047 – 0.080)	0.094 (0.071 – 0.118)
	Merl-R18	0.007 (0.005 – 0.011)	0.931 (0.912 – 0.947)	0.436 (0.372 – 0.505)	0.984 (0.983 – 0.984)	0.014 (0.011 – 0.017)	0.054 (0.045 – 0.064)
2° 1 Type AVB	Random Init-R18	0.005 (0.003 – 0.007)	0.774 (0.748 – 0.800)	0.097 (0.066 – 0.135)	0.998 (0.997 – 0.998)	0.026 (0.018 – 0.037)	0.038 (0.025 – 0.052)
	ECGCLIP-R18	0.007 (0.005 – 0.010)	0.766 (0.739 – 0.795)	0.181 (0.142 – 0.227)	0.996 (0.995 – 0.996)	0.029 (0.022 – 0.036)	0.052 (0.040 – 0.065)
	ECGCLIP-R34	0.007 (0.005 – 0.010)	0.741 (0.708 – 0.772)	0.171 (0.130 – 0.217)	0.997 (0.997 – 0.997)	0.036 (0.027 – 0.046)	0.057 (0.043 – 0.073)
	Merl-R18	0.007 (0.004 – 0.010)	0.801 (0.777 – 0.825)	0.040 (0.019 – 0.065)	0.999 (0.999 – 0.999)	0.030 (0.014 – 0.047)	0.031 (0.014 – 0.048)
2° 2 Type AVB	Random Init-R18	0.005 (0.004 – 0.007)	0.610 (0.595 – 0.625)	0.012 (0.007 – 0.018)	0.999 (0.999 – 0.999)	0.018 (0.011 – 0.026)	0.020 (0.012 – 0.029)
	ECGCLIP-R18	0.006 (0.004 – 0.008)	0.611 (0.597 – 0.625)	0.016 (0.011 – 0.023)	0.999 (0.999 – 0.999)	0.024 (0.015 – 0.033)	0.026 (0.016 – 0.037)
	ECGCLIP-R34	0.006 (0.004 – 0.010)	0.501 (0.485 – 0.515)	0.014 (0.009 – 0.020)	0.999 (0.999 – 0.999)	0.021 (0.013 – 0.029)	0.022 (0.013 – 0.032)
	Merl-R18	0.004 (0.003 – 0.006)	0.503 (0.489 – 0.518)	0.007 (0.003 – 0.011)	1.000 (1.000 – 1.000)	0.012 (0.005 – 0.020)	0.019 (0.008 – 0.032)
LVH	Random Init-R18	0.401 (0.397 – 0.405)	0.800 (0.798 – 0.802)	0.170 (0.167 – 0.173)	0.992 (0.992 – 0.992)	0.272 (0.268 – 0.276)	0.313 (0.309 – 0.317)
	ECGCLIP-R18	0.416 (0.412 – 0.419)	0.820 (0.818 – 0.821)	0.162 (0.159 – 0.164)	0.993 (0.992 – 0.993)	0.262 (0.258 – 0.266)	0.307 (0.303 – 0.311)
	ECGCLIP-R34	0.403 (0.399 – 0.407)	0.816 (0.815 – 0.818)	0.185 (0.182 – 0.187)	0.990 (0.990 – 0.990)	0.288 (0.284 – 0.291)	0.316 (0.312 – 0.320)
	Merl-R18	0.256 (0.253 – 0.259)	0.711 (0.709 – 0.713)	0.153 (0.150 – 0.155)	0.976 (0.975 – 0.976)	0.220 (0.216 – 0.223)	0.200 (0.196 – 0.203)
LAH	Random Init-R18	0.085 (0.083 – 0.088)	0.832 (0.829 – 0.834)	0.334 (0.327 – 0.341)	0.931 (0.930 – 0.931)	0.151 (0.148 – 0.155)	0.147 (0.144 – 0.151)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RVH	ECGCLIP-R18	0.085 (0.083 - 0.087)	0.836 (0.833 - 0.838)	0.291 (0.284 - 0.298)	0.940 (0.939 - 0.941)	0.147 (0.143 - 0.151)	0.138 (0.133 - 0.142)
	ECGCLIP-R34	0.082 (0.080 - 0.084)	0.831 (0.828 - 0.833)	0.382 (0.375 - 0.389)	0.913 (0.913 - 0.914)	0.146 (0.143 - 0.149)	0.149 (0.145 - 0.153)
	Merl-R18	0.181 (0.176 - 0.187)	0.870 (0.867 - 0.872)	0.353 (0.346 - 0.360)	0.970 (0.969 - 0.970)	0.262 (0.256 - 0.267)	0.250 (0.244 - 0.255)
	Random Init-R18	0.067 (0.063 - 0.071)	0.804 (0.799 - 0.809)	0.112 (0.105 - 0.119)	0.991 (0.991 - 0.992)	0.117 (0.110 - 0.124)	0.108 (0.101 - 0.115)
	ECGCLIP-R18	0.080 (0.076 - 0.085)	0.832 (0.827 - 0.836)	0.085 (0.079 - 0.092)	0.996 (0.996 - 0.996)	0.114 (0.106 - 0.122)	0.114 (0.106 - 0.123)
	ECGCLIP-R34	0.076 (0.072 - 0.080)	0.823 (0.818 - 0.828)	0.093 (0.087 - 0.100)	0.995 (0.995 - 0.995)	0.118 (0.110 - 0.126)	0.115 (0.107 - 0.124)
	Merl-R18	0.102 (0.097 - 0.108)	0.838 (0.833 - 0.842)	0.061 (0.056 - 0.067)	0.999 (0.998 - 0.999)	0.102 (0.093 - 0.111)	0.133 (0.122 - 0.143)
	Random Init-R18	0.002 (0.002 - 0.002)	0.813 (0.791 - 0.834)	0.012 (0.003 - 0.026)	0.999 (0.999 - 0.999)	0.006 (0.001 - 0.013)	0.006 (0.001 - 0.014)
RAH	ECGCLIP-R18	0.005 (0.004 - 0.006)	0.886 (0.870 - 0.901)	0.033 (0.016 - 0.054)	0.999 (0.999 - 0.999)	0.021 (0.010 - 0.035)	0.022 (0.010 - 0.037)
	ECGCLIP-R34	0.006 (0.004 - 0.007)	0.886 (0.865 - 0.905)	0.006 (0.000 - 0.015)	1.000 (1.000 - 1.000)	0.007 (0.000 - 0.019)	0.007 (0.000 - 0.019)
	Merl-R18	0.008 (0.006 - 0.010)	0.917 (0.902 - 0.929)	0.012 (0.003 - 0.026)	1.000 (1.000 - 1.000)	0.013 (0.003 - 0.028)	0.013 (0.003 - 0.027)
	Random Init-R18	0.006 (0.004 - 0.009)	0.911 (0.893 - 0.927)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.000 - 0.001)	0.000 (0.000 - 0.000)
DC	ECGCLIP-R18	0.002 (0.002 - 0.004)	0.814 (0.784 - 0.845)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.002 (0.002 - 0.004)	0.802 (0.767 - 0.838)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	Merl-R18	0.001 (0.001 - 0.002)	0.758 (0.725 - 0.791)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.000 - 0.001)	0.000 (0.000 - 0.000)
BVH	Random Init-R18	0.093 (0.090 - 0.097)	0.796 (0.793 - 0.799)	0.001 (0.001 - 0.002)	1.000 (1.000 - 1.000)	0.002 (0.001 - 0.003)	0.021 (0.014 - 0.029)
	ECGCLIP-R18	0.088 (0.086 - 0.091)	0.790 (0.787 - 0.793)	0.001 (0.001 - 0.002)	1.000 (1.000 - 1.000)	0.003 (0.002 - 0.004)	0.018 (0.012 - 0.026)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
VAT	ECGCLIP-R34	0.084 (0.082 - 0.087)	0.796 (0.794 - 0.800)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.001)	0.004 (-0.001 - 0.010)
	Merl-R18	0.059 (0.057 - 0.060)	0.721 (0.718 - 0.725)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.001)	0.004 (-0.001 - 0.009)
	Random Init-R18	0.018 (0.015 - 0.022)	0.986 (0.983 - 0.989)	0.490 (0.425 - 0.553)	0.993 (0.993 - 0.994)	0.043 (0.036 - 0.051)	0.104 (0.089 - 0.119)
	ECGCLIP-R18	0.015 (0.012 - 0.017)	0.985 (0.982 - 0.988)	0.406 (0.348 - 0.469)	0.993 (0.993 - 0.994)	0.036 (0.030 - 0.044)	0.087 (0.073 - 0.102)
	ECGCLIP-R34	0.011 (0.009 - 0.012)	0.984 (0.981 - 0.986)	0.172 (0.126 - 0.222)	0.994 (0.994 - 0.994)	0.017 (0.012 - 0.022)	0.038 (0.026 - 0.049)
	Merl-R18	0.008 (0.006 - 0.010)	0.956 (0.946 - 0.965)	0.577 (0.514 - 0.642)	0.978 (0.978 - 0.978)	0.016 (0.014 - 0.019)	0.067 (0.058 - 0.075)
DDD	Random Init-R18	0.082 (0.067 - 0.101)	0.993 (0.990 - 0.996)	0.750 (0.698 - 0.799)	0.996 (0.995 - 0.996)	0.113 (0.100 - 0.127)	0.213 (0.195 - 0.231)
	ECGCLIP-R18	0.089 (0.074 - 0.108)	0.996 (0.995 - 0.997)	0.851 (0.808 - 0.892)	0.993 (0.993 - 0.993)	0.083 (0.074 - 0.093)	0.192 (0.177 - 0.207)
	ECGCLIP-R34	0.100 (0.083 - 0.122)	0.996 (0.995 - 0.997)	0.851 (0.811 - 0.892)	0.992 (0.992 - 0.992)	0.076 (0.068 - 0.085)	0.184 (0.170 - 0.197)
	Merl-R18	0.094 (0.076 - 0.116)	0.995 (0.992 - 0.996)	0.778 (0.731 - 0.823)	0.994 (0.993 - 0.994)	0.082 (0.073 - 0.093)	0.183 (0.168 - 0.198)
HK	Random Init-R18	0.010 (0.009 - 0.011)	0.762 (0.750 - 0.774)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.004 (0.004 - 0.005)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.010 (0.009 - 0.011)	0.725 (0.710 - 0.739)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.004 (0.004 - 0.005)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.008 (0.007 - 0.008)	0.716 (0.702 - 0.730)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.004 (0.004 - 0.005)	0.000 (0.000 - 0.000)
	Merl-R18	0.006 (0.005 - 0.006)	0.696 (0.683 - 0.708)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.004 (0.004 - 0.005)	0.000 (0.000 - 0.000)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S21: Detailed task-specific performance comparison on the external UKB cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STEMI	Random Init-R18	0.009 (0.006 - 0.021)	0.640 (0.602 - 0.677)	0.005 (0.000 - 0.015)	1.000 (1.000 - 1.000)	0.009 (0.000 - 0.029)	0.047 (0.000 - 0.116)
	ECGCLIP-R18	0.017 (0.008 - 0.033)	0.682 (0.645 - 0.717)	0.005 (0.000 - 0.015)	1.000 (1.000 - 1.000)	0.009 (0.000 - 0.029)	0.033 (-0.001 - 0.096)
	ECGCLIP-R34	0.020 (0.012 - 0.033)	0.732 (0.698 - 0.766)	0.005 (0.000 - 0.015)	1.000 (1.000 - 1.000)	0.009 (0.000 - 0.029)	0.039 (-0.001 - 0.102)
	Merl-R18	0.015 (0.008 - 0.026)	0.690 (0.654 - 0.726)	0.005 (0.000 - 0.015)	1.000 (1.000 - 1.000)	0.009 (0.000 - 0.029)	0.033 (-0.001 - 0.097)
NSR	Random Init-R18	0.975 (0.973 - 0.976)	0.655 (0.647 - 0.663)	0.543 (0.539 - 0.547)	0.720 (0.704 - 0.738)	0.697 (0.694 - 0.701)	0.115 (0.107 - 0.123)
	ECGCLIP-R18	0.975 (0.973 - 0.976)	0.652 (0.644 - 0.661)	0.537 (0.533 - 0.541)	0.718 (0.701 - 0.735)	0.692 (0.688 - 0.696)	0.111 (0.104 - 0.119)
	ECGCLIP-R34	0.978 (0.977 - 0.979)	0.692 (0.684 - 0.700)	0.546 (0.542 - 0.551)	0.836 (0.822 - 0.850)	0.703 (0.699 - 0.706)	0.167 (0.159 - 0.174)
	Merl-R18	0.973 (0.972 - 0.975)	0.640 (0.631 - 0.648)	0.549 (0.544 - 0.553)	0.675 (0.657 - 0.693)	0.701 (0.697 - 0.704)	0.098 (0.090 - 0.106)
SBrad	Random Init-R18	0.886 (0.882 - 0.889)	0.905 (0.902 - 0.908)	0.603 (0.597 - 0.610)	0.944 (0.941 - 0.946)	0.721 (0.716 - 0.726)	0.593 (0.587 - 0.600)
	ECGCLIP-R18	0.887 (0.883 - 0.891)	0.910 (0.907 - 0.912)	0.632 (0.625 - 0.638)	0.937 (0.934 - 0.939)	0.739 (0.734 - 0.744)	0.607 (0.600 - 0.614)
	ECGCLIP-R34	0.984 (0.983 - 0.985)	0.986 (0.986 - 0.987)	0.908 (0.904 - 0.912)	0.967 (0.965 - 0.969)	0.932 (0.930 - 0.934)	0.880 (0.876 - 0.884)
	Merl-R18	0.842 (0.837 - 0.847)	0.872 (0.869 - 0.875)	0.496 (0.489 - 0.502)	0.945 (0.942 - 0.948)	0.634 (0.628 - 0.640)	0.505 (0.498 - 0.512)
VPB	Random Init-R18	0.346 (0.319 - 0.373)	0.829 (0.817 - 0.841)	0.256 (0.234 - 0.279)	0.997 (0.996 - 0.997)	0.372 (0.344 - 0.400)	0.408 (0.380 - 0.434)
	ECGCLIP-R18	0.445 (0.417 - 0.476)	0.881 (0.871 - 0.890)	0.322 (0.298 - 0.348)	0.997 (0.997 - 0.998)	0.454 (0.427 - 0.482)	0.490 (0.463 - 0.515)
	ECGCLIP-R34	0.616 (0.588 - 0.645)	0.940 (0.932 - 0.947)	0.538 (0.510 - 0.566)	0.995 (0.995 - 0.996)	0.635 (0.611 - 0.659)	0.637 (0.614 - 0.660)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
APB	Merl-R18	0.367 (0.339 - 0.394)	0.858 (0.847 - 0.869)	0.245 (0.222 - 0.267)	0.997 (0.997 - 0.998)	0.366 (0.337 - 0.393)	0.413 (0.385 - 0.438)
	Random Init-R18	0.029 (0.027 - 0.032)	0.613 (0.598 - 0.630)	0.522 (0.492 - 0.551)	0.629 (0.624 - 0.633)	0.055 (0.050 - 0.059)	0.045 (0.035 - 0.053)
	ECGCLIP-R18	0.034 (0.031 - 0.038)	0.644 (0.628 - 0.660)	0.483 (0.453 - 0.513)	0.711 (0.707 - 0.715)	0.064 (0.059 - 0.069)	0.061 (0.051 - 0.070)
	ECGCLIP-R34	0.087 (0.074 - 0.103)	0.743 (0.730 - 0.757)	0.469 (0.439 - 0.500)	0.836 (0.833 - 0.839)	0.102 (0.094 - 0.111)	0.116 (0.104 - 0.128)
SArr	Merl-R18	0.032 (0.029 - 0.035)	0.642 (0.625 - 0.657)	0.513 (0.483 - 0.547)	0.682 (0.678 - 0.686)	0.062 (0.057 - 0.067)	0.059 (0.050 - 0.070)
	Random Init-R18	0.046 (0.044 - 0.049)	0.495 (0.485 - 0.507)	0.001 (0.000 - 0.002)	0.998 (0.997 - 0.998)	0.002 (0.000 - 0.004)	-0.006 (-0.011 - 0.000)
	ECGCLIP-R18	0.046 (0.043 - 0.048)	0.496 (0.486 - 0.508)	0.003 (0.001 - 0.005)	0.996 (0.996 - 0.997)	0.005 (0.002 - 0.009)	-0.004 (-0.011 - 0.003)
	ECGCLIP-R34	0.045 (0.043 - 0.048)	0.507 (0.496 - 0.518)	0.069 (0.060 - 0.080)	0.905 (0.903 - 0.908)	0.046 (0.040 - 0.053)	-0.019 (-0.025 - -0.011)
AF	Merl-R18	0.048 (0.045 - 0.051)	0.512 (0.501 - 0.523)	0.003 (0.001 - 0.005)	0.996 (0.996 - 0.997)	0.005 (0.002 - 0.009)	-0.004 (-0.010 - 0.004)
	Random Init-R18	0.853 (0.829 - 0.880)	0.984 (0.979 - 0.989)	0.602 (0.567 - 0.640)	0.999 (0.999 - 1.000)	0.737 (0.709 - 0.766)	0.753 (0.729 - 0.779)
	ECGCLIP-R18	0.865 (0.843 - 0.888)	0.985 (0.981 - 0.990)	0.644 (0.610 - 0.676)	1.000 (0.999 - 1.000)	0.769 (0.742 - 0.793)	0.781 (0.758 - 0.803)
	ECGCLIP-R34	0.795 (0.767 - 0.821)	0.976 (0.969 - 0.982)	0.627 (0.592 - 0.660)	0.999 (0.999 - 0.999)	0.743 (0.716 - 0.769)	0.753 (0.728 - 0.777)
AFL	Merl-R18	0.840 (0.816 - 0.868)	0.983 (0.978 - 0.988)	0.627 (0.592 - 0.662)	0.999 (0.999 - 1.000)	0.750 (0.722 - 0.777)	0.763 (0.738 - 0.787)
	Random Init-R18	0.455 (0.323 - 0.586)	0.910 (0.865 - 0.948)	0.354 (0.237 - 0.477)	1.000 (1.000 - 1.000)	0.479 (0.342 - 0.600)	0.512 (0.384 - 0.626)
	ECGCLIP-R18	0.452 (0.316 - 0.580)	0.916 (0.871 - 0.953)	0.400 (0.281 - 0.519)	1.000 (1.000 - 1.000)	0.525 (0.396 - 0.636)	0.553 (0.431 - 0.657)
	ECGCLIP-R34	0.449 (0.316 - 0.586)	0.899 (0.845 - 0.943)	0.415 (0.289 - 0.533)	1.000 (1.000 - 1.000)	0.540 (0.405 - 0.651)	0.566 (0.445 - 0.667)
	Merl-R18	0.429 (0.291 - 0.560)	0.879 (0.823 - 0.931)	0.385 (0.265 - 0.500)	1.000 (1.000 - 1.000)	0.505 (0.372 - 0.618)	0.531 (0.408 - 0.640)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
JPB	Random Init-R18	0.011 (0.010 - 0.013)	0.553 (0.526 - 0.579)	0.130 (0.101 - 0.161)	0.926 (0.924 - 0.929)	0.028 (0.021 - 0.035)	0.020 (0.010 - 0.031)
	ECGCLIP-R18	0.012 (0.010 - 0.014)	0.574 (0.548 - 0.601)	0.029 (0.014 - 0.045)	0.986 (0.985 - 0.987)	0.023 (0.012 - 0.036)	0.012 (0.001 - 0.026)
	ECGCLIP-R34	0.018 (0.015 - 0.022)	0.634 (0.607 - 0.660)	0.064 (0.042 - 0.088)	0.985 (0.984 - 0.986)	0.047 (0.031 - 0.065)	0.038 (0.021 - 0.056)
	Merl-R18	0.013 (0.011 - 0.016)	0.588 (0.561 - 0.613)	0.022 (0.009 - 0.037)	0.992 (0.991 - 0.993)	0.023 (0.010 - 0.039)	0.015 (0.002 - 0.030)
AJR	Random Init-R18	0.025 (0.018 - 0.035)	0.848 (0.821 - 0.875)	0.277 (0.214 - 0.345)	0.969 (0.968 - 0.971)	0.057 (0.043 - 0.073)	0.085 (0.063 - 0.108)
	ECGCLIP-R18	0.026 (0.019 - 0.034)	0.860 (0.834 - 0.886)	0.272 (0.210 - 0.335)	0.976 (0.975 - 0.978)	0.070 (0.054 - 0.089)	0.097 (0.072 - 0.123)
	ECGCLIP-R34	0.030 (0.022 - 0.046)	0.880 (0.857 - 0.902)	0.212 (0.158 - 0.273)	0.982 (0.981 - 0.984)	0.070 (0.052 - 0.093)	0.087 (0.063 - 0.116)
	Merl-R18	0.025 (0.019 - 0.033)	0.849 (0.821 - 0.876)	0.174 (0.123 - 0.229)	0.982 (0.981 - 0.983)	0.056 (0.039 - 0.075)	0.069 (0.046 - 0.093)
JTach	Random Init-R18	0.026 (0.020 - 0.035)	0.883 (0.862 - 0.903)	0.011 (0.000 - 0.027)	0.999 (0.999 - 1.000)	0.018 (0.000 - 0.045)	0.023 (-0.002 - 0.062)
	ECGCLIP-R18	0.030 (0.022 - 0.041)	0.889 (0.868 - 0.908)	0.016 (0.000 - 0.035)	1.000 (0.999 - 1.000)	0.029 (0.000 - 0.061)	0.042 (-0.001 - 0.091)
	ECGCLIP-R34	0.033 (0.025 - 0.045)	0.894 (0.870 - 0.915)	0.016 (0.000 - 0.035)	1.000 (0.999 - 1.000)	0.029 (0.000 - 0.061)	0.043 (-0.001 - 0.093)
	Merl-R18	0.026 (0.020 - 0.037)	0.886 (0.867 - 0.903)	0.022 (0.005 - 0.044)	0.999 (0.999 - 1.000)	0.037 (0.009 - 0.074)	0.050 (0.010 - 0.098)
JEB	Random Init-R18	0.025 (0.019 - 0.035)	0.836 (0.808 - 0.867)	0.109 (0.066 - 0.157)	0.993 (0.992 - 0.994)	0.073 (0.044 - 0.105)	0.072 (0.042 - 0.105)
	ECGCLIP-R18	0.025 (0.019 - 0.034)	0.837 (0.810 - 0.866)	0.016 (0.000 - 0.036)	0.998 (0.998 - 0.998)	0.021 (0.000 - 0.046)	0.019 (-0.003 - 0.045)
	ECGCLIP-R34	0.035 (0.025 - 0.049)	0.866 (0.842 - 0.892)	0.027 (0.006 - 0.053)	0.998 (0.998 - 0.999)	0.038 (0.008 - 0.072)	0.039 (0.006 - 0.075)
	Merl-R18	0.024 (0.018 - 0.032)	0.834 (0.805 - 0.865)	0.033 (0.011 - 0.060)	0.997 (0.996 - 0.997)	0.034 (0.011 - 0.062)	0.030 (0.008 - 0.059)
RBBB	Random Init-R18	0.793 (0.780 - 0.806)	0.974 (0.971 - 0.977)	0.541 (0.521 - 0.561)	0.996 (0.995 - 0.996)	0.667 (0.649 - 0.684)	0.674 (0.656 - 0.689)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
1° AVB	ECGCLIP-R18	0.785 (0.771 – 0.798)	0.974 (0.971 – 0.976)	0.578 (0.558 – 0.597)	0.992 (0.992 – 0.993)	0.671 (0.655 – 0.687)	0.667 (0.650 – 0.682)
	ECGCLIP-R34	0.793 (0.780 – 0.805)	0.975 (0.972 – 0.978)	0.608 (0.588 – 0.627)	0.991 (0.990 – 0.992)	0.684 (0.668 – 0.699)	0.675 (0.659 – 0.690)
	Merl-R18	0.806 (0.793 – 0.819)	0.978 (0.975 – 0.980)	0.560 (0.539 – 0.579)	0.995 (0.994 – 0.995)	0.675 (0.658 – 0.691)	0.677 (0.661 – 0.692)
	Random Init-R18	0.856 (0.844 – 0.867)	0.988 (0.987 – 0.989)	0.771 (0.756 – 0.786)	0.986 (0.985 – 0.987)	0.773 (0.760 – 0.785)	0.759 (0.746 – 0.771)
	ECGCLIP-R18	0.879 (0.869 – 0.889)	0.990 (0.989 – 0.991)	0.856 (0.843 – 0.868)	0.981 (0.980 – 0.982)	0.791 (0.780 – 0.801)	0.779 (0.768 – 0.790)
	ECGCLIP-R34	0.906 (0.896 – 0.914)	0.992 (0.991 – 0.993)	0.864 (0.851 – 0.876)	0.985 (0.984 – 0.986)	0.823 (0.813 – 0.833)	0.812 (0.802 – 0.823)
	Merl-R18	0.868 (0.856 – 0.879)	0.989 (0.988 – 0.990)	0.844 (0.830 – 0.857)	0.980 (0.979 – 0.981)	0.780 (0.769 – 0.791)	0.768 (0.756 – 0.779)
	Random Init-R18	0.096 (0.076 – 0.118)	0.847 (0.831 – 0.862)	0.008 (0.002 – 0.016)	1.000 (1.000 – 1.000)	0.015 (0.004 – 0.031)	0.055 (0.013 – 0.101)
	ECGCLIP-R18	0.111 (0.087 – 0.137)	0.857 (0.842 – 0.871)	0.018 (0.008 – 0.030)	1.000 (1.000 – 1.000)	0.034 (0.015 – 0.058)	0.105 (0.053 – 0.155)
	ECGCLIP-R34	0.137 (0.108 – 0.168)	0.889 (0.877 – 0.900)	0.045 (0.028 – 0.064)	1.000 (0.999 – 1.000)	0.084 (0.053 – 0.116)	0.159 (0.109 – 0.205)
QT Prolong	Merl-R18	0.060 (0.046 – 0.080)	0.797 (0.779 – 0.815)	0.004 (0.000 – 0.010)	1.000 (1.000 – 1.000)	0.008 (0.000 – 0.019)	0.039 (-0.001 – 0.085)
	Random Init-R18	0.240 (0.201 – 0.284)	0.931 (0.918 – 0.944)	0.000 (0.000 – 0.000)	1.000 (1.000 – 1.000)	0.000 (0.000 – 0.000)	0.000 (0.000 – 0.000)
	ECGCLIP-R18	0.268 (0.227 – 0.315)	0.931 (0.919 – 0.944)	0.000 (0.000 – 0.000)	1.000 (1.000 – 1.000)	0.000 (0.000 – 0.000)	0.000 (0.000 – 0.000)
	ECGCLIP-R34	0.322 (0.278 – 0.370)	0.945 (0.935 – 0.955)	0.005 (0.000 – 0.012)	1.000 (1.000 – 1.000)	0.009 (0.000 – 0.023)	0.068 (0.000 – 0.108)
ER	Merl-R18	0.230 (0.192 – 0.273)	0.917 (0.904 – 0.931)	0.000 (0.000 – 0.000)	1.000 (1.000 – 1.000)	0.000 (0.000 – 0.000)	0.000 (0.000 – 0.000)
	Random Init-R18	0.496 (0.457 – 0.532)	0.974 (0.970 – 0.978)	0.070 (0.051 – 0.089)	1.000 (1.000 – 1.000)	0.129 (0.096 – 0.161)	0.237 (0.192 – 0.274)
	ECGCLIP-R18	0.591 (0.555 – 0.628)	0.983 (0.979 – 0.986)	0.099 (0.079 – 0.121)	1.000 (1.000 – 1.000)	0.179 (0.145 – 0.215)	0.299 (0.261 – 0.334)
LAFB							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LBBB	ECGCLIP-R34	0.625 (0.586 - 0.660)	0.985 (0.981 - 0.988)	0.148 (0.124 - 0.174)	1.000 (0.999 - 1.000)	0.251 (0.214 - 0.289)	0.346 (0.308 - 0.382)
	Merl-R18	0.564 (0.527 - 0.601)	0.981 (0.977 - 0.984)	0.123 (0.102 - 0.148)	1.000 (0.999 - 1.000)	0.215 (0.181 - 0.252)	0.320 (0.284 - 0.358)
	Random Init-R18	0.879 (0.847 - 0.909)	0.984 (0.977 - 0.991)	0.246 (0.208 - 0.283)	1.000 (1.000 - 1.000)	0.390 (0.340 - 0.436)	0.482 (0.441 - 0.518)
	ECGCLIP-R18	0.890 (0.860 - 0.919)	0.986 (0.979 - 0.992)	0.329 (0.289 - 0.368)	1.000 (1.000 - 1.000)	0.490 (0.444 - 0.532)	0.560 (0.521 - 0.595)
	ECGCLIP-R34	0.887 (0.857 - 0.917)	0.981 (0.972 - 0.988)	0.485 (0.442 - 0.526)	1.000 (1.000 - 1.000)	0.644 (0.603 - 0.680)	0.679 (0.643 - 0.710)
	Merl-R18	0.895 (0.867 - 0.922)	0.982 (0.973 - 0.990)	0.248 (0.212 - 0.283)	1.000 (1.000 - 1.000)	0.394 (0.347 - 0.439)	0.486 (0.445 - 0.522)
IVB	Random Init-R18	0.057 (0.045 - 0.076)	0.838 (0.813 - 0.863)	0.036 (0.017 - 0.058)	0.998 (0.998 - 0.999)	0.054 (0.026 - 0.086)	0.059 (0.027 - 0.096)
	ECGCLIP-R18	0.052 (0.041 - 0.068)	0.815 (0.790 - 0.842)	0.039 (0.019 - 0.062)	0.998 (0.997 - 0.998)	0.056 (0.028 - 0.088)	0.059 (0.027 - 0.094)
	ECGCLIP-R34	0.050 (0.039 - 0.064)	0.812 (0.783 - 0.838)	0.036 (0.016 - 0.057)	0.997 (0.997 - 0.998)	0.049 (0.022 - 0.079)	0.049 (0.020 - 0.081)
	Merl-R18	0.053 (0.043 - 0.067)	0.851 (0.828 - 0.874)	0.032 (0.014 - 0.054)	0.998 (0.998 - 0.998)	0.048 (0.020 - 0.078)	0.051 (0.020 - 0.083)
Short PR	Random Init-R18	0.349 (0.301 - 0.400)	0.938 (0.924 - 0.950)	0.034 (0.018 - 0.053)	1.000 (1.000 - 1.000)	0.066 (0.035 - 0.100)	0.164 (0.106 - 0.216)
	ECGCLIP-R18	0.346 (0.298 - 0.396)	0.933 (0.918 - 0.947)	0.027 (0.014 - 0.044)	1.000 (1.000 - 1.000)	0.053 (0.028 - 0.084)	0.153 (0.100 - 0.202)
	ECGCLIP-R34	0.359 (0.310 - 0.410)	0.937 (0.924 - 0.950)	0.094 (0.068 - 0.122)	1.000 (1.000 - 1.000)	0.167 (0.123 - 0.212)	0.263 (0.209 - 0.315)
	Merl-R18	0.331 (0.281 - 0.378)	0.936 (0.922 - 0.949)	0.043 (0.024 - 0.063)	1.000 (1.000 - 1.000)	0.082 (0.047 - 0.117)	0.177 (0.118 - 0.227)
LVH	Random Init-R18	0.088 (0.075 - 0.103)	0.715 (0.697 - 0.732)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.105 (0.090 - 0.123)	0.732 (0.715 - 0.748)	0.001 (0.000 - 0.003)	1.000 (1.000 - 1.000)	0.002 (0.000 - 0.006)	0.031 (0.000 - 0.055)
	ECGCLIP-R34	0.085 (0.074 - 0.099)	0.741 (0.726 - 0.756)	0.001 (0.000 - 0.003)	1.000 (1.000 - 1.000)	0.002 (0.000 - 0.006)	0.031 (0.000 - 0.055)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAH	Merl-R18	0.076 (0.066 - 0.088)	0.729 (0.712 - 0.745)	0.002 (0.000 - 0.005)	1.000 (1.000 - 1.000)	0.004 (0.000 - 0.010)	0.044 (0.000 - 0.069)
	Random Init-R18	0.397 (0.380 - 0.415)	0.883 (0.877 - 0.888)	0.002 (0.001 - 0.003)	1.000 (1.000 - 1.000)	0.004 (0.001 - 0.007)	0.033 (0.012 - 0.051)
	ECGCLIP-R18	0.620 (0.603 - 0.636)	0.942 (0.938 - 0.945)	0.001 (0.000 - 0.003)	1.000 (1.000 - 1.000)	0.002 (0.001 - 0.005)	0.030 (0.008 - 0.048)
	ECGCLIP-R34	0.515 (0.497 - 0.531)	0.925 (0.921 - 0.929)	0.006 (0.003 - 0.009)	1.000 (1.000 - 1.000)	0.012 (0.007 - 0.017)	0.068 (0.048 - 0.086)
RVH	Merl-R18	0.416 (0.399 - 0.432)	0.890 (0.885 - 0.895)	0.003 (0.001 - 0.005)	1.000 (1.000 - 1.000)	0.005 (0.002 - 0.009)	0.045 (0.025 - 0.064)
	Random Init-R18	0.011 (0.007 - 0.018)	0.799 (0.751 - 0.847)	0.000 (0.000 - 0.000)	1.000 (0.999 - 1.000)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)
	ECGCLIP-R18	0.016 (0.009 - 0.026)	0.855 (0.818 - 0.889)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)
	ECGCLIP-R34	0.012 (0.007 - 0.020)	0.841 (0.803 - 0.876)	0.000 (0.000 - 0.000)	1.000 (0.999 - 1.000)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)
RAH	Merl-R18	0.012 (0.007 - 0.018)	0.826 (0.786 - 0.866)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
	Random Init-R18	0.016 (0.008 - 0.032)	0.737 (0.688 - 0.786)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.027 (0.015 - 0.047)	0.854 (0.820 - 0.883)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.018 (0.011 - 0.030)	0.842 (0.806 - 0.874)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Merl-R18	0.023 (0.012 - 0.044)	0.783 (0.736 - 0.824)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S22: Detailed task-specific performance comparison on the external Chapman cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STEMI	Random Init-R18	0.087 (0.055 – 0.138)	0.894 (0.852 – 0.930)	0.130 (0.074 – 0.187)	0.997 (0.997 – 0.998)	0.119 (0.067 – 0.169)	0.117 (0.065 – 0.168)
	ECGCLIP-R18	0.095 (0.065 – 0.132)	0.925 (0.890 – 0.954)	0.163 (0.099 – 0.225)	0.997 (0.996 – 0.997)	0.141 (0.085 – 0.193)	0.140 (0.083 – 0.192)
	ECGCLIP-R34	0.116 (0.083 – 0.155)	0.943 (0.912 – 0.968)	0.341 (0.258 – 0.421)	0.994 (0.994 – 0.995)	0.200 (0.146 – 0.250)	0.217 (0.160 – 0.270)
	Merl-R18	0.081 (0.050 – 0.122)	0.915 (0.879 – 0.945)	0.179 (0.113 – 0.240)	0.996 (0.996 – 0.997)	0.140 (0.087 – 0.189)	0.141 (0.087 – 0.189)
SBrad	Random Init-R18	0.977 (0.973 – 0.981)	0.993 (0.992 – 0.994)	0.913 (0.909 – 0.917)	0.986 (0.984 – 0.987)	0.942 (0.940 – 0.945)	0.912 (0.908 – 0.916)
	ECGCLIP-R18	0.990 (0.988 – 0.991)	0.996 (0.995 – 0.996)	0.912 (0.908 – 0.916)	0.990 (0.988 – 0.991)	0.945 (0.943 – 0.948)	0.917 (0.913 – 0.920)
	ECGCLIP-R34	0.996 (0.995 – 0.997)	0.998 (0.998 – 0.998)	0.954 (0.951 – 0.957)	0.993 (0.992 – 0.994)	0.971 (0.969 – 0.972)	0.955 (0.952 – 0.957)
	Merl-R18	0.983 (0.980 – 0.986)	0.994 (0.993 – 0.994)	0.864 (0.858 – 0.869)	0.990 (0.989 – 0.991)	0.918 (0.915 – 0.921)	0.879 (0.875 – 0.884)
STach	Random Init-R18	0.971 (0.968 – 0.975)	0.995 (0.994 – 0.995)	0.859 (0.851 – 0.867)	0.993 (0.993 – 0.994)	0.907 (0.902 – 0.913)	0.893 (0.887 – 0.899)
	ECGCLIP-R18	0.980 (0.977 – 0.982)	0.996 (0.995 – 0.996)	0.843 (0.834 – 0.851)	0.994 (0.993 – 0.995)	0.900 (0.895 – 0.906)	0.885 (0.879 – 0.891)
	ECGCLIP-R34	0.975 (0.973 – 0.977)	0.994 (0.994 – 0.995)	0.842 (0.833 – 0.850)	0.997 (0.996 – 0.997)	0.906 (0.900 – 0.911)	0.893 (0.887 – 0.898)
	Merl-R18	0.973 (0.970 – 0.976)	0.994 (0.994 – 0.995)	0.832 (0.823 – 0.840)	0.994 (0.993 – 0.994)	0.892 (0.886 – 0.898)	0.876 (0.870 – 0.883)
VPB	Random Init-R18	0.132 (0.108 – 0.161)	0.916 (0.896 – 0.934)	0.641 (0.585 – 0.693)	0.966 (0.964 – 0.967)	0.194 (0.171 – 0.217)	0.260 (0.234 – 0.286)
	ECGCLIP-R18	0.161 (0.132 – 0.191)	0.922 (0.899 – 0.941)	0.686 (0.633 – 0.737)	0.968 (0.966 – 0.970)	0.217 (0.193 – 0.241)	0.288 (0.260 – 0.314)
	ECGCLIP-R34	0.210 (0.175 – 0.243)	0.935 (0.915 – 0.952)	0.754 (0.708 – 0.803)	0.969 (0.968 – 0.971)	0.244 (0.216 – 0.269)	0.322 (0.295 – 0.347)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
APB	Merl-R18	0.135 (0.110 - 0.162)	0.917 (0.897 - 0.935)	0.631 (0.573 - 0.680)	0.969 (0.967 - 0.970)	0.205 (0.180 - 0.229)	0.268 (0.240 - 0.294)
	Random Init-R18	0.000 (0.000 - 0.001)	0.876 (0.813 - 0.953)	0.667 (0.000 - 1.000)	0.828 (0.824 - 0.831)	0.001 (0.000 - 0.001)	0.011 (-0.003 - 0.023)
	ECGCLIP-R18	0.000 (0.000 - 0.001)	0.891 (0.822 - 0.960)	0.667 (0.000 - 1.000)	0.882 (0.880 - 0.885)	0.001 (0.000 - 0.002)	0.014 (-0.002 - 0.029)
	ECGCLIP-R34	0.001 (0.000 - 0.002)	0.958 (0.942 - 0.970)	1.000 (1.000 - 1.000)	0.943 (0.941 - 0.945)	0.002 (0.001 - 0.005)	0.033 (0.019 - 0.050)
SArr	Merl-R18	0.000 (0.000 - 0.001)	0.911 (0.825 - 0.965)	0.667 (0.000 - 1.000)	0.850 (0.846 - 0.853)	0.001 (0.000 - 0.001)	0.012 (-0.003 - 0.025)
	Random Init-R18	0.231 (0.216 - 0.247)	0.789 (0.780 - 0.797)	0.230 (0.214 - 0.246)	0.974 (0.973 - 0.976)	0.277 (0.260 - 0.294)	0.249 (0.231 - 0.267)
	ECGCLIP-R18	0.234 (0.218 - 0.250)	0.793 (0.785 - 0.802)	0.277 (0.260 - 0.294)	0.966 (0.964 - 0.967)	0.299 (0.282 - 0.315)	0.262 (0.244 - 0.278)
	ECGCLIP-R34	0.258 (0.241 - 0.274)	0.831 (0.824 - 0.837)	0.376 (0.357 - 0.394)	0.944 (0.942 - 0.946)	0.324 (0.309 - 0.339)	0.281 (0.266 - 0.297)
AF	Merl-R18	0.235 (0.220 - 0.252)	0.800 (0.792 - 0.807)	0.338 (0.320 - 0.357)	0.945 (0.943 - 0.947)	0.299 (0.283 - 0.314)	0.254 (0.237 - 0.270)
	Random Init-R18	0.204 (0.192 - 0.217)	0.906 (0.902 - 0.911)	0.655 (0.632 - 0.678)	0.895 (0.892 - 0.898)	0.312 (0.297 - 0.326)	0.322 (0.308 - 0.336)
	ECGCLIP-R18	0.201 (0.189 - 0.215)	0.904 (0.900 - 0.910)	0.625 (0.601 - 0.649)	0.901 (0.898 - 0.903)	0.309 (0.294 - 0.324)	0.315 (0.299 - 0.330)
	ECGCLIP-R34	0.192 (0.180 - 0.205)	0.892 (0.886 - 0.898)	0.614 (0.590 - 0.639)	0.904 (0.901 - 0.906)	0.310 (0.296 - 0.325)	0.314 (0.298 - 0.330)
AFL	Merl-R18	0.203 (0.192 - 0.217)	0.905 (0.900 - 0.910)	0.684 (0.663 - 0.707)	0.890 (0.887 - 0.893)	0.315 (0.301 - 0.329)	0.330 (0.317 - 0.345)
	Random Init-R18	0.641 (0.631 - 0.652)	0.909 (0.906 - 0.911)	0.128 (0.121 - 0.135)	0.996 (0.995 - 0.996)	0.223 (0.212 - 0.234)	0.294 (0.283 - 0.305)
	ECGCLIP-R18	0.654 (0.644 - 0.664)	0.907 (0.905 - 0.910)	0.143 (0.135 - 0.150)	0.996 (0.996 - 0.997)	0.246 (0.234 - 0.257)	0.319 (0.308 - 0.329)
	ECGCLIP-R34	0.664 (0.654 - 0.674)	0.905 (0.902 - 0.908)	0.143 (0.135 - 0.150)	0.997 (0.997 - 0.998)	0.247 (0.235 - 0.258)	0.327 (0.316 - 0.336)
	Merl-R18	0.632 (0.622 - 0.643)	0.905 (0.902 - 0.908)	0.126 (0.120 - 0.133)	0.996 (0.996 - 0.997)	0.221 (0.210 - 0.232)	0.298 (0.287 - 0.309)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
JPB	Random Init-R18	0.000 (0.000 - 0.001)	0.532 (0.397 - 0.676)	0.000 (0.000 - 0.000)	0.941 (0.939 - 0.943)	0.000 (0.000 - 0.000)	-0.004 (-0.005 - -0.003)
	ECGCLIP-R18	0.000 (0.000 - 0.001)	0.626 (0.465 - 0.767)	0.000 (0.000 - 0.000)	0.973 (0.972 - 0.975)	0.000 (0.000 - 0.000)	-0.003 (-0.003 - -0.002)
	ECGCLIP-R34	0.001 (0.000 - 0.002)	0.689 (0.531 - 0.833)	0.091 (0.000 - 0.286)	0.986 (0.985 - 0.987)	0.003 (0.000 - 0.010)	0.010 (-0.002 - 0.036)
	Merl-R18	0.000 (0.000 - 0.001)	0.564 (0.431 - 0.703)	0.000 (0.000 - 0.000)	0.984 (0.983 - 0.985)	0.000 (0.000 - 0.000)	-0.002 (-0.003 - -0.001)
ATach	Random Init-R18	0.087 (0.070 - 0.107)	0.922 (0.908 - 0.935)	0.578 (0.524 - 0.634)	0.954 (0.952 - 0.956)	0.137 (0.119 - 0.154)	0.199 (0.176 - 0.222)
	ECGCLIP-R18	0.137 (0.110 - 0.170)	0.938 (0.926 - 0.949)	0.524 (0.469 - 0.582)	0.976 (0.975 - 0.977)	0.204 (0.177 - 0.233)	0.248 (0.218 - 0.280)
	ECGCLIP-R34	0.195 (0.160 - 0.238)	0.958 (0.949 - 0.965)	0.571 (0.520 - 0.623)	0.983 (0.981 - 0.984)	0.271 (0.240 - 0.302)	0.311 (0.279 - 0.342)
	Merl-R18	0.095 (0.078 - 0.116)	0.936 (0.927 - 0.945)	0.527 (0.474 - 0.586)	0.966 (0.965 - 0.968)	0.160 (0.137 - 0.182)	0.212 (0.186 - 0.238)
JEB	Random Init-R18	0.034 (0.023 - 0.048)	0.931 (0.901 - 0.955)	0.373 (0.266 - 0.481)	0.989 (0.988 - 0.990)	0.092 (0.061 - 0.125)	0.137 (0.093 - 0.181)
	ECGCLIP-R18	0.038 (0.026 - 0.055)	0.940 (0.915 - 0.959)	0.200 (0.115 - 0.292)	0.993 (0.993 - 0.994)	0.078 (0.045 - 0.115)	0.096 (0.055 - 0.139)
	ECGCLIP-R34	0.051 (0.032 - 0.093)	0.958 (0.939 - 0.972)	0.120 (0.051 - 0.191)	0.995 (0.994 - 0.996)	0.057 (0.025 - 0.091)	0.064 (0.026 - 0.103)
	Merl-R18	0.035 (0.023 - 0.051)	0.926 (0.895 - 0.952)	0.160 (0.082 - 0.240)	0.995 (0.994 - 0.995)	0.075 (0.039 - 0.113)	0.085 (0.043 - 0.130)
SVT	Random Init-R18	0.759 (0.726 - 0.791)	0.973 (0.964 - 0.980)	0.494 (0.457 - 0.529)	0.999 (0.999 - 0.999)	0.636 (0.602 - 0.666)	0.660 (0.628 - 0.686)
	ECGCLIP-R18	0.777 (0.745 - 0.809)	0.978 (0.970 - 0.984)	0.498 (0.461 - 0.532)	0.999 (0.999 - 0.999)	0.644 (0.611 - 0.674)	0.669 (0.641 - 0.696)
	ECGCLIP-R34	0.821 (0.791 - 0.849)	0.986 (0.980 - 0.991)	0.502 (0.467 - 0.537)	0.999 (0.999 - 1.000)	0.653 (0.621 - 0.683)	0.681 (0.653 - 0.707)
	Merl-R18	0.753 (0.718 - 0.785)	0.980 (0.974 - 0.984)	0.470 (0.434 - 0.505)	0.999 (0.999 - 0.999)	0.615 (0.581 - 0.645)	0.643 (0.613 - 0.670)
VEB	Random Init-R18	0.028 (0.016 - 0.050)	0.897 (0.850 - 0.937)	0.175 (0.089 - 0.267)	0.996 (0.996 - 0.997)	0.091 (0.044 - 0.142)	0.101 (0.050 - 0.157)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RBBB	ECGCLIP-R18	0.031 (0.019 - 0.052)	0.918 (0.877 - 0.953)	0.063 (0.014 - 0.125)	0.999 (0.998 - 0.999)	0.066 (0.015 - 0.128)	0.065 (0.014 - 0.127)
	ECGCLIP-R34	0.025 (0.017 - 0.039)	0.934 (0.907 - 0.957)	0.032 (0.000 - 0.079)	0.998 (0.997 - 0.998)	0.024 (0.000 - 0.060)	0.023 (-0.002 - 0.060)
	Merl-R18	0.029 (0.018 - 0.049)	0.902 (0.854 - 0.943)	0.254 (0.153 - 0.364)	0.992 (0.991 - 0.992)	0.070 (0.039 - 0.103)	0.099 (0.056 - 0.142)
	Random Init-R18	0.350 (0.316 - 0.383)	0.982 (0.978 - 0.985)	0.949 (0.932 - 0.964)	0.959 (0.957 - 0.960)	0.397 (0.374 - 0.420)	0.477 (0.458 - 0.495)
	ECGCLIP-R18	0.336 (0.306 - 0.370)	0.982 (0.978 - 0.985)	0.940 (0.921 - 0.957)	0.960 (0.959 - 0.962)	0.404 (0.382 - 0.427)	0.481 (0.462 - 0.500)
	ECGCLIP-R34	0.338 (0.305 - 0.371)	0.982 (0.978 - 0.985)	0.946 (0.929 - 0.963)	0.958 (0.956 - 0.960)	0.394 (0.372 - 0.416)	0.474 (0.455 - 0.492)
	Merl-R18	0.348 (0.315 - 0.383)	0.982 (0.978 - 0.985)	0.948 (0.931 - 0.964)	0.960 (0.958 - 0.962)	0.404 (0.382 - 0.425)	0.482 (0.464 - 0.500)
	Random Init-R18	0.652 (0.628 - 0.677)	0.888 (0.875 - 0.901)	0.573 (0.546 - 0.599)	0.993 (0.992 - 0.993)	0.634 (0.611 - 0.656)	0.627 (0.605 - 0.650)
	ECGCLIP-R18	0.659 (0.635 - 0.684)	0.880 (0.867 - 0.894)	0.607 (0.579 - 0.634)	0.991 (0.990 - 0.992)	0.645 (0.624 - 0.667)	0.636 (0.614 - 0.658)
	ECGCLIP-R34	0.668 (0.645 - 0.692)	0.878 (0.864 - 0.892)	0.572 (0.544 - 0.599)	0.994 (0.993 - 0.995)	0.650 (0.628 - 0.671)	0.647 (0.624 - 0.668)
1° AVB	Merl-R18	0.653 (0.627 - 0.679)	0.892 (0.880 - 0.905)	0.628 (0.601 - 0.653)	0.989 (0.988 - 0.990)	0.638 (0.617 - 0.660)	0.627 (0.605 - 0.649)
	Random Init-R18	0.116 (0.095 - 0.145)	0.914 (0.902 - 0.926)	0.218 (0.181 - 0.258)	0.991 (0.990 - 0.991)	0.190 (0.159 - 0.224)	0.184 (0.153 - 0.217)
	ECGCLIP-R18	0.152 (0.124 - 0.189)	0.929 (0.919 - 0.938)	0.254 (0.211 - 0.296)	0.991 (0.990 - 0.992)	0.221 (0.185 - 0.257)	0.215 (0.180 - 0.252)
	ECGCLIP-R34	0.199 (0.162 - 0.243)	0.941 (0.932 - 0.949)	0.294 (0.253 - 0.344)	0.990 (0.989 - 0.991)	0.244 (0.210 - 0.283)	0.240 (0.205 - 0.279)
QT Prolong	Merl-R18	0.132 (0.106 - 0.165)	0.913 (0.902 - 0.925)	0.170 (0.135 - 0.211)	0.993 (0.992 - 0.994)	0.175 (0.140 - 0.214)	0.168 (0.133 - 0.207)
	Random Init-R18	0.241 (0.202 - 0.285)	0.951 (0.940 - 0.962)	0.266 (0.222 - 0.310)	0.996 (0.996 - 0.997)	0.308 (0.261 - 0.352)	0.308 (0.261 - 0.350)
	ECGCLIP-R18	0.268 (0.224 - 0.318)	0.958 (0.949 - 0.967)	0.192 (0.152 - 0.233)	0.998 (0.998 - 0.999)	0.275 (0.224 - 0.324)	0.302 (0.249 - 0.351)
ER							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAFB	ECGCLIP-R34	0.272 (0.225 - 0.318)	0.963 (0.954 - 0.970)	0.321 (0.272 - 0.368)	0.995 (0.994 - 0.996)	0.329 (0.282 - 0.370)	0.323 (0.277 - 0.365)
	Merl-R18	0.218 (0.177 - 0.264)	0.934 (0.921 - 0.946)	0.219 (0.178 - 0.262)	0.997 (0.996 - 0.997)	0.271 (0.225 - 0.315)	0.275 (0.227 - 0.320)
	Random Init-R18	0.071 (0.049 - 0.102)	0.828 (0.801 - 0.854)	0.208 (0.157 - 0.264)	0.992 (0.991 - 0.993)	0.151 (0.113 - 0.192)	0.151 (0.112 - 0.193)
	ECGCLIP-R18	0.070 (0.050 - 0.099)	0.840 (0.813 - 0.866)	0.167 (0.122 - 0.215)	0.994 (0.993 - 0.995)	0.144 (0.106 - 0.185)	0.140 (0.102 - 0.181)
	ECGCLIP-R34	0.071 (0.052 - 0.096)	0.847 (0.819 - 0.872)	0.200 (0.151 - 0.251)	0.993 (0.992 - 0.993)	0.155 (0.118 - 0.194)	0.153 (0.117 - 0.194)
	Merl-R18	0.065 (0.046 - 0.093)	0.809 (0.776 - 0.839)	0.158 (0.114 - 0.207)	0.994 (0.994 - 0.995)	0.144 (0.105 - 0.183)	0.139 (0.101 - 0.180)
IVB	Random Init-R18	0.029 (0.025 - 0.035)	0.574 (0.553 - 0.597)	0.032 (0.021 - 0.046)	0.995 (0.994 - 0.995)	0.049 (0.032 - 0.069)	0.047 (0.028 - 0.069)
	ECGCLIP-R18	0.031 (0.027 - 0.037)	0.589 (0.568 - 0.612)	0.027 (0.017 - 0.040)	0.995 (0.995 - 0.996)	0.042 (0.026 - 0.061)	0.042 (0.023 - 0.065)
	ECGCLIP-R34	0.033 (0.029 - 0.039)	0.612 (0.591 - 0.635)	0.029 (0.018 - 0.041)	0.996 (0.995 - 0.996)	0.045 (0.028 - 0.063)	0.046 (0.026 - 0.068)
	Merl-R18	0.029 (0.025 - 0.034)	0.574 (0.550 - 0.596)	0.026 (0.016 - 0.039)	0.995 (0.995 - 0.996)	0.040 (0.024 - 0.058)	0.039 (0.020 - 0.060)
VPE	Random Init-R18	0.625 (0.505 - 0.752)	0.984 (0.971 - 0.994)	0.514 (0.403 - 0.635)	1.000 (1.000 - 1.000)	0.627 (0.529 - 0.727)	0.642 (0.552 - 0.739)
	ECGCLIP-R18	0.640 (0.523 - 0.756)	0.985 (0.974 - 0.994)	0.542 (0.426 - 0.661)	1.000 (1.000 - 1.000)	0.645 (0.537 - 0.739)	0.656 (0.554 - 0.747)
	ECGCLIP-R34	0.652 (0.534 - 0.769)	0.990 (0.982 - 0.996)	0.486 (0.375 - 0.605)	1.000 (1.000 - 1.000)	0.609 (0.495 - 0.709)	0.629 (0.523 - 0.720)
	Merl-R18	0.610 (0.498 - 0.720)	0.986 (0.977 - 0.994)	0.500 (0.389 - 0.615)	1.000 (1.000 - 1.000)	0.605 (0.504 - 0.704)	0.618 (0.519 - 0.714)
3° AVB	Random Init-R18	0.033 (0.025 - 0.043)	0.684 (0.646 - 0.717)	0.128 (0.092 - 0.164)	0.989 (0.988 - 0.990)	0.095 (0.068 - 0.123)	0.090 (0.062 - 0.119)
	ECGCLIP-R18	0.049 (0.032 - 0.070)	0.696 (0.656 - 0.727)	0.150 (0.114 - 0.188)	0.987 (0.986 - 0.988)	0.102 (0.076 - 0.129)	0.099 (0.072 - 0.127)
	ECGCLIP-R34	0.067 (0.045 - 0.094)	0.746 (0.716 - 0.773)	0.081 (0.052 - 0.111)	0.997 (0.996 - 0.997)	0.104 (0.068 - 0.138)	0.103 (0.066 - 0.139)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
2° 1 Type AVB	Merl-R18	0.028 (0.022 - 0.035)	0.706 (0.670 - 0.737)	0.138 (0.102 - 0.176)	0.985 (0.984 - 0.986)	0.084 (0.062 - 0.108)	0.081 (0.057 - 0.108)
	Random Init-R18	0.035 (0.025 - 0.049)	0.707 (0.675 - 0.735)	0.021 (0.006 - 0.035)	0.999 (0.999 - 1.000)	0.037 (0.011 - 0.063)	0.060 (0.017 - 0.101)
	ECGCLIP-R18	0.037 (0.026 - 0.052)	0.664 (0.628 - 0.699)	0.029 (0.013 - 0.048)	0.999 (0.999 - 0.999)	0.050 (0.022 - 0.080)	0.066 (0.029 - 0.107)
	ECGCLIP-R34	0.038 (0.024 - 0.057)	0.680 (0.649 - 0.710)	0.044 (0.023 - 0.069)	0.999 (0.999 - 0.999)	0.076 (0.041 - 0.116)	0.106 (0.059 - 0.158)
2° 2 Type AVB	Merl-R18	0.030 (0.022 - 0.041)	0.673 (0.639 - 0.703)	0.009 (0.000 - 0.020)	0.999 (0.999 - 1.000)	0.016 (0.000 - 0.036)	0.030 (-0.002 - 0.068)
	Random Init-R18	0.013 (0.010 - 0.019)	0.591 (0.558 - 0.624)	0.016 (0.003 - 0.031)	0.997 (0.997 - 0.998)	0.023 (0.005 - 0.044)	0.021 (0.001 - 0.044)
	ECGCLIP-R18	0.019 (0.013 - 0.030)	0.620 (0.587 - 0.654)	0.013 (0.003 - 0.027)	0.999 (0.999 - 1.000)	0.024 (0.006 - 0.050)	0.041 (0.009 - 0.082)
	ECGCLIP-R34	0.017 (0.010 - 0.032)	0.443 (0.408 - 0.483)	0.016 (0.003 - 0.031)	1.000 (1.000 - 1.000)	0.031 (0.006 - 0.058)	0.064 (0.014 - 0.117)
LVH	Merl-R18	0.011 (0.008 - 0.016)	0.494 (0.459 - 0.533)	0.003 (0.000 - 0.011)	0.999 (0.998 - 0.999)	0.005 (0.000 - 0.018)	0.004 (-0.003 - 0.021)
	Random Init-R18	0.358 (0.320 - 0.395)	0.971 (0.967 - 0.974)	0.930 (0.911 - 0.949)	0.904 (0.901 - 0.907)	0.218 (0.204 - 0.232)	0.319 (0.306 - 0.333)
	ECGCLIP-R18	0.336 (0.299 - 0.374)	0.970 (0.966 - 0.973)	0.936 (0.918 - 0.954)	0.904 (0.901 - 0.906)	0.219 (0.204 - 0.233)	0.321 (0.307 - 0.335)
	ECGCLIP-R34	0.359 (0.322 - 0.396)	0.972 (0.969 - 0.975)	0.952 (0.936 - 0.967)	0.888 (0.885 - 0.891)	0.198 (0.184 - 0.211)	0.303 (0.290 - 0.315)
RVH	Merl-R18	0.363 (0.325 - 0.402)	0.971 (0.967 - 0.975)	0.943 (0.925 - 0.959)	0.890 (0.887 - 0.893)	0.198 (0.185 - 0.211)	0.302 (0.289 - 0.315)
	Random Init-R18	0.260 (0.189 - 0.349)	0.984 (0.973 - 0.992)	0.700 (0.612 - 0.780)	0.993 (0.992 - 0.993)	0.296 (0.241 - 0.352)	0.360 (0.305 - 0.415)
	ECGCLIP-R18	0.264 (0.191 - 0.350)	0.984 (0.975 - 0.991)	0.618 (0.521 - 0.709)	0.995 (0.994 - 0.995)	0.323 (0.267 - 0.378)	0.365 (0.306 - 0.421)
	ECGCLIP-R34	0.285 (0.209 - 0.382)	0.986 (0.979 - 0.992)	0.591 (0.500 - 0.680)	0.995 (0.994 - 0.995)	0.318 (0.259 - 0.377)	0.356 (0.295 - 0.415)
	Merl-R18	0.261 (0.188 - 0.348)	0.981 (0.966 - 0.991)	0.536 (0.439 - 0.624)	0.996 (0.996 - 0.997)	0.354 (0.287 - 0.420)	0.375 (0.305 - 0.439)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RAH	Random Init-R18	0.251 (0.121 - 0.425)	0.981 (0.968 - 0.992)	0.583 (0.410 - 0.743)	0.997 (0.997 - 0.998)	0.241 (0.156 - 0.324)	0.297 (0.203 - 0.380)
	ECGCLIP-R18	0.251 (0.127 - 0.410)	0.977 (0.950 - 0.995)	0.556 (0.382 - 0.710)	0.998 (0.997 - 0.998)	0.260 (0.167 - 0.346)	0.306 (0.204 - 0.394)
	ECGCLIP-R34	0.313 (0.161 - 0.495)	0.957 (0.895 - 0.994)	0.472 (0.305 - 0.639)	0.999 (0.999 - 0.999)	0.358 (0.227 - 0.471)	0.368 (0.241 - 0.484)
	Merl-R18	0.253 (0.122 - 0.423)	0.972 (0.949 - 0.991)	0.417 (0.258 - 0.586)	0.999 (0.999 - 0.999)	0.357 (0.222 - 0.488)	0.360 (0.227 - 0.490)
LK	Random Init-R18	0.018 (0.013 - 0.024)	0.849 (0.814 - 0.878)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.006 (0.005 - 0.007)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.019 (0.014 - 0.026)	0.856 (0.821 - 0.884)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.006 (0.005 - 0.007)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.025 (0.017 - 0.036)	0.878 (0.850 - 0.902)	0.044 (0.014 - 0.079)	0.998 (0.998 - 0.999)	0.057 (0.019 - 0.098)	0.057 (0.017 - 0.102)
	Merl-R18	0.018 (0.013 - 0.025)	0.836 (0.801 - 0.868)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.006 (0.005 - 0.007)	0.000 (0.000 - 0.000)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S23: Detailed task-specific performance comparison on the external PTB-XL cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STEMI	Random Init-R18	0.755 (0.744 – 0.765)	0.871 (0.865 – 0.876)	0.036 (0.031 – 0.041)	1.000 (1.000 – 1.000)	0.070 (0.060 – 0.079)	0.167 (0.154 – 0.177)
	ECGCLIP-R18	0.759 (0.749 – 0.770)	0.872 (0.866 – 0.877)	0.039 (0.034 – 0.044)	1.000 (1.000 – 1.000)	0.076 (0.066 – 0.085)	0.173 (0.161 – 0.185)
	ECGCLIP-R34	0.762 (0.751 – 0.772)	0.882 (0.877 – 0.888)	0.052 (0.046 – 0.058)	1.000 (1.000 – 1.000)	0.099 (0.089 – 0.110)	0.198 (0.186 – 0.210)
	Merl-R18	0.779 (0.769 – 0.788)	0.883 (0.878 – 0.888)	0.065 (0.058 – 0.071)	1.000 (1.000 – 1.000)	0.122 (0.109 – 0.133)	0.222 (0.210 – 0.233)
NSR	Random Init-R18	0.551 (0.541 – 0.562)	0.615 (0.608 – 0.622)	0.783 (0.775 – 0.792)	0.336 (0.327 – 0.344)	0.593 (0.586 – 0.600)	0.131 (0.118 – 0.144)
	ECGCLIP-R18	0.533 (0.522 – 0.544)	0.610 (0.603 – 0.618)	0.758 (0.750 – 0.767)	0.357 (0.349 – 0.366)	0.586 (0.579 – 0.594)	0.125 (0.112 – 0.138)
	ECGCLIP-R34	0.513 (0.503 – 0.524)	0.594 (0.586 – 0.601)	0.757 (0.749 – 0.766)	0.352 (0.343 – 0.360)	0.583 (0.577 – 0.591)	0.117 (0.104 – 0.130)
	Merl-R18	0.552 (0.541 – 0.562)	0.616 (0.609 – 0.624)	0.774 (0.766 – 0.782)	0.346 (0.337 – 0.354)	0.591 (0.584 – 0.598)	0.130 (0.118 – 0.143)
SBrad	Random Init-R18	0.427 (0.386 – 0.469)	0.957 (0.951 – 0.962)	0.929 (0.909 – 0.949)	0.870 (0.866 – 0.874)	0.298 (0.280 – 0.316)	0.374 (0.358 – 0.389)
	ECGCLIP-R18	0.342 (0.306 – 0.377)	0.949 (0.943 – 0.955)	0.934 (0.915 – 0.953)	0.864 (0.860 – 0.869)	0.290 (0.274 – 0.308)	0.368 (0.352 – 0.384)
	ECGCLIP-R34	0.291 (0.261 – 0.323)	0.946 (0.940 – 0.952)	0.948 (0.930 – 0.965)	0.855 (0.850 – 0.859)	0.280 (0.264 – 0.298)	0.361 (0.347 – 0.376)
	Merl-R18	0.366 (0.330 – 0.404)	0.951 (0.945 – 0.957)	0.915 (0.893 – 0.935)	0.875 (0.870 – 0.879)	0.302 (0.284 – 0.320)	0.375 (0.359 – 0.391)
STach	Random Init-R18	0.874 (0.848 – 0.897)	0.991 (0.987 – 0.994)	0.835 (0.811 – 0.860)	0.995 (0.994 – 0.996)	0.847 (0.829 – 0.864)	0.841 (0.822 – 0.860)
	ECGCLIP-R18	0.880 (0.856 – 0.904)	0.990 (0.986 – 0.994)	0.856 (0.833 – 0.880)	0.995 (0.994 – 0.996)	0.861 (0.844 – 0.878)	0.855 (0.838 – 0.873)
	ECGCLIP-R34	0.895 (0.869 – 0.918)	0.992 (0.987 – 0.995)	0.898 (0.878 – 0.919)	0.994 (0.993 – 0.995)	0.876 (0.861 – 0.892)	0.871 (0.856 – 0.888)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
VPB	Merl-R18	0.875 (0.849 – 0.899)	0.989 (0.984 – 0.993)	0.856 (0.832 – 0.878)	0.994 (0.993 – 0.995)	0.851 (0.834 – 0.868)	0.845 (0.827 – 0.863)
	Random Init-R18	0.712 (0.679 – 0.742)	0.968 (0.964 – 0.973)	0.763 (0.737 – 0.789)	0.978 (0.976 – 0.980)	0.707 (0.687 – 0.727)	0.692 (0.671 – 0.713)
	ECGCLIP-R18	0.761 (0.730 – 0.789)	0.975 (0.971 – 0.979)	0.794 (0.771 – 0.818)	0.984 (0.982 – 0.985)	0.760 (0.741 – 0.779)	0.746 (0.727 – 0.766)
	ECGCLIP-R34	0.835 (0.811 – 0.858)	0.986 (0.983 – 0.989)	0.899 (0.883 – 0.916)	0.984 (0.983 – 0.986)	0.823 (0.807 – 0.839)	0.816 (0.800 – 0.832)
APB	Merl-R18	0.722 (0.690 – 0.751)	0.968 (0.963 – 0.973)	0.747 (0.723 – 0.771)	0.983 (0.982 – 0.985)	0.729 (0.708 – 0.749)	0.714 (0.692 – 0.736)
	Random Init-R18	0.067 (0.056 – 0.082)	0.793 (0.773 – 0.813)	0.525 (0.477 – 0.572)	0.852 (0.847 – 0.856)	0.111 (0.097 – 0.124)	0.139 (0.120 – 0.157)
	ECGCLIP-R18	0.093 (0.078 – 0.114)	0.833 (0.813 – 0.850)	0.477 (0.428 – 0.527)	0.901 (0.897 – 0.905)	0.140 (0.123 – 0.159)	0.165 (0.142 – 0.187)
	ECGCLIP-R34	0.197 (0.165 – 0.236)	0.897 (0.882 – 0.910)	0.558 (0.506 – 0.604)	0.949 (0.946 – 0.952)	0.260 (0.233 – 0.286)	0.286 (0.257 – 0.313)
SArr	Merl-R18	0.074 (0.062 – 0.090)	0.809 (0.788 – 0.829)	0.475 (0.427 – 0.525)	0.886 (0.882 – 0.890)	0.125 (0.109 – 0.141)	0.149 (0.128 – 0.170)
	Random Init-R18	0.068 (0.059 – 0.077)	0.638 (0.617 – 0.660)	0.060 (0.045 – 0.078)	0.983 (0.981 – 0.984)	0.078 (0.058 – 0.100)	0.057 (0.037 – 0.081)
	ECGCLIP-R18	0.068 (0.060 – 0.078)	0.642 (0.621 – 0.662)	0.093 (0.072 – 0.115)	0.972 (0.970 – 0.974)	0.101 (0.078 – 0.124)	0.071 (0.048 – 0.094)
	ECGCLIP-R34	0.070 (0.062 – 0.080)	0.657 (0.636 – 0.676)	0.175 (0.147 – 0.201)	0.940 (0.936 – 0.943)	0.124 (0.104 – 0.142)	0.086 (0.064 – 0.106)
AF	Merl-R18	0.064 (0.057 – 0.074)	0.626 (0.604 – 0.647)	0.079 (0.060 – 0.099)	0.973 (0.970 – 0.975)	0.087 (0.065 – 0.108)	0.057 (0.035 – 0.079)
	Random Init-R18	0.758 (0.734 – 0.784)	0.969 (0.964 – 0.974)	0.637 (0.612 – 0.662)	0.990 (0.989 – 0.992)	0.721 (0.701 – 0.741)	0.711 (0.690 – 0.730)
	ECGCLIP-R18	0.772 (0.747 – 0.797)	0.964 (0.958 – 0.970)	0.575 (0.549 – 0.601)	0.993 (0.992 – 0.994)	0.688 (0.666 – 0.709)	0.684 (0.663 – 0.704)
	ECGCLIP-R34	0.776 (0.751 – 0.802)	0.966 (0.960 – 0.971)	0.591 (0.567 – 0.617)	0.993 (0.992 – 0.994)	0.703 (0.682 – 0.723)	0.699 (0.679 – 0.719)
	Merl-R18	0.774 (0.750 – 0.798)	0.965 (0.959 – 0.971)	0.635 (0.610 – 0.660)	0.991 (0.990 – 0.992)	0.723 (0.704 – 0.742)	0.714 (0.694 – 0.732)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AFL	Random Init-R18	0.622 (0.500 - 0.733)	0.987 (0.971 - 0.997)	0.781 (0.679 - 0.872)	0.996 (0.996 - 0.997)	0.548 (0.462 - 0.626)	0.572 (0.489 - 0.647)
	ECGCLIP-R18	0.641 (0.519 - 0.755)	0.990 (0.977 - 0.998)	0.836 (0.743 - 0.914)	0.996 (0.995 - 0.997)	0.540 (0.458 - 0.615)	0.575 (0.504 - 0.642)
	ECGCLIP-R34	0.628 (0.507 - 0.738)	0.992 (0.986 - 0.998)	0.863 (0.778 - 0.938)	0.996 (0.995 - 0.997)	0.560 (0.479 - 0.633)	0.596 (0.527 - 0.662)
	Merl-R18	0.589 (0.465 - 0.708)	0.990 (0.977 - 0.998)	0.822 (0.727 - 0.905)	0.996 (0.995 - 0.997)	0.538 (0.456 - 0.617)	0.572 (0.494 - 0.644)
AJR	Random Init-R18	0.002 (0.001 - 0.003)	0.712 (0.648 - 0.775)	0.000 (0.000 - 0.000)	0.994 (0.993 - 0.995)	0.000 (0.000 - 0.000)	-0.003 (-0.003 - -0.002)
	ECGCLIP-R18	0.005 (0.002 - 0.012)	0.794 (0.736 - 0.847)	0.111 (0.000 - 0.250)	0.993 (0.992 - 0.994)	0.032 (0.000 - 0.073)	0.043 (-0.003 - 0.101)
	ECGCLIP-R34	0.008 (0.004 - 0.018)	0.887 (0.843 - 0.923)	0.074 (0.000 - 0.192)	0.995 (0.994 - 0.996)	0.029 (0.000 - 0.074)	0.034 (-0.003 - 0.089)
	Merl-R18	0.003 (0.001 - 0.004)	0.762 (0.701 - 0.815)	0.000 (0.000 - 0.000)	0.996 (0.995 - 0.997)	0.000 (0.000 - 0.000)	-0.002 (-0.003 - -0.002)
ATach	Random Init-R18	0.117 (0.032 - 0.255)	0.940 (0.887 - 0.982)	0.593 (0.400 - 0.786)	0.986 (0.984 - 0.987)	0.091 (0.051 - 0.134)	0.168 (0.107 - 0.230)
	ECGCLIP-R18	0.065 (0.011 - 0.167)	0.860 (0.758 - 0.939)	0.296 (0.133 - 0.476)	0.991 (0.990 - 0.992)	0.071 (0.027 - 0.117)	0.106 (0.042 - 0.171)
	ECGCLIP-R34	0.027 (0.013 - 0.053)	0.955 (0.924 - 0.978)	0.296 (0.136 - 0.500)	0.990 (0.989 - 0.992)	0.066 (0.026 - 0.117)	0.102 (0.042 - 0.172)
	Merl-R18	0.051 (0.016 - 0.164)	0.927 (0.854 - 0.979)	0.333 (0.158 - 0.531)	0.991 (0.990 - 0.992)	0.080 (0.035 - 0.132)	0.121 (0.053 - 0.192)
JTach	Random Init-R18	0.024 (0.011 - 0.053)	0.936 (0.897 - 0.970)	0.111 (0.000 - 0.250)	0.997 (0.997 - 0.998)	0.068 (0.000 - 0.146)	0.072 (-0.002 - 0.159)
	ECGCLIP-R18	0.043 (0.015 - 0.108)	0.948 (0.908 - 0.976)	0.185 (0.048 - 0.345)	0.997 (0.996 - 0.998)	0.105 (0.024 - 0.195)	0.115 (0.026 - 0.216)
	ECGCLIP-R34	0.026 (0.009 - 0.069)	0.923 (0.877 - 0.962)	0.148 (0.033 - 0.296)	0.997 (0.996 - 0.997)	0.077 (0.018 - 0.159)	0.086 (0.018 - 0.175)
	Merl-R18	0.024 (0.010 - 0.058)	0.942 (0.904 - 0.969)	0.111 (0.000 - 0.250)	0.997 (0.996 - 0.998)	0.061 (0.000 - 0.130)	0.067 (-0.002 - 0.145)
SVT	Random Init-R18	0.057 (0.027 - 0.094)	0.643 (0.596 - 0.688)	0.028 (0.006 - 0.054)	1.000 (0.999 - 1.000)	0.053 (0.011 - 0.100)	0.111 (0.029 - 0.188)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RBBB	ECGCLIP-R18	0.062 (0.032 - 0.102)	0.715 (0.674 - 0.754)	0.034 (0.011 - 0.064)	1.000 (1.000 - 1.000)	0.063 (0.021 - 0.118)	0.140 (0.054 - 0.222)
	ECGCLIP-R34	0.064 (0.033 - 0.104)	0.705 (0.661 - 0.749)	0.039 (0.012 - 0.070)	1.000 (1.000 - 1.000)	0.074 (0.024 - 0.129)	0.164 (0.072 - 0.243)
	Merl-R18	0.066 (0.035 - 0.106)	0.678 (0.633 - 0.722)	0.039 (0.012 - 0.070)	1.000 (1.000 - 1.000)	0.073 (0.023 - 0.127)	0.150 (0.062 - 0.226)
	Random Init-R18	0.880 (0.863 - 0.895)	0.979 (0.975 - 0.983)	0.650 (0.626 - 0.673)	0.996 (0.995 - 0.996)	0.763 (0.745 - 0.779)	0.761 (0.743 - 0.776)
	ECGCLIP-R18	0.882 (0.866 - 0.897)	0.979 (0.975 - 0.982)	0.623 (0.597 - 0.646)	0.996 (0.996 - 0.997)	0.747 (0.727 - 0.766)	0.748 (0.729 - 0.765)
	ECGCLIP-R34	0.884 (0.868 - 0.899)	0.978 (0.974 - 0.982)	0.640 (0.615 - 0.661)	0.996 (0.995 - 0.997)	0.759 (0.741 - 0.776)	0.759 (0.741 - 0.774)
	Merl-R18	0.879 (0.863 - 0.894)	0.978 (0.973 - 0.981)	0.642 (0.617 - 0.664)	0.996 (0.995 - 0.996)	0.757 (0.739 - 0.775)	0.755 (0.738 - 0.772)
	Random Init-R18	0.507 (0.471 - 0.548)	0.955 (0.948 - 0.962)	0.677 (0.643 - 0.711)	0.968 (0.966 - 0.970)	0.536 (0.508 - 0.562)	0.527 (0.499 - 0.552)
	ECGCLIP-R18	0.507 (0.470 - 0.545)	0.958 (0.952 - 0.964)	0.704 (0.672 - 0.736)	0.966 (0.964 - 0.968)	0.541 (0.513 - 0.566)	0.535 (0.508 - 0.560)
	ECGCLIP-R34	0.511 (0.476 - 0.547)	0.963 (0.958 - 0.968)	0.681 (0.647 - 0.713)	0.971 (0.969 - 0.973)	0.557 (0.532 - 0.583)	0.547 (0.521 - 0.573)
1° AVB	Merl-R18	0.518 (0.481 - 0.555)	0.960 (0.954 - 0.966)	0.721 (0.691 - 0.754)	0.965 (0.963 - 0.968)	0.545 (0.520 - 0.570)	0.542 (0.517 - 0.566)
	Random Init-R18	0.079 (0.045 - 0.132)	0.862 (0.835 - 0.893)	0.179 (0.119 - 0.257)	0.992 (0.991 - 0.994)	0.139 (0.092 - 0.195)	0.137 (0.089 - 0.195)
	ECGCLIP-R18	0.096 (0.056 - 0.150)	0.884 (0.858 - 0.910)	0.248 (0.179 - 0.330)	0.993 (0.991 - 0.994)	0.190 (0.136 - 0.249)	0.190 (0.134 - 0.251)
QT Prolong	ECGCLIP-R34	0.101 (0.064 - 0.160)	0.886 (0.859 - 0.909)	0.256 (0.184 - 0.337)	0.991 (0.989 - 0.992)	0.171 (0.122 - 0.220)	0.175 (0.125 - 0.227)
	Merl-R18	0.084 (0.050 - 0.136)	0.871 (0.842 - 0.899)	0.214 (0.149 - 0.293)	0.993 (0.992 - 0.994)	0.170 (0.118 - 0.227)	0.168 (0.115 - 0.227)
	Random Init-R18	0.696 (0.671 - 0.717)	0.953 (0.948 - 0.957)	0.273 (0.251 - 0.293)	0.996 (0.996 - 0.997)	0.414 (0.387 - 0.438)	0.465 (0.440 - 0.487)
LAFB	ECGCLIP-R18	0.679 (0.656 - 0.703)	0.941 (0.935 - 0.946)	0.176 (0.159 - 0.194)	0.998 (0.998 - 0.999)	0.295 (0.269 - 0.319)	0.382 (0.358 - 0.404)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LBBB	ECGCLIP-R34	0.706 (0.684 - 0.729)	0.951 (0.945 - 0.956)	0.227 (0.207 - 0.246)	0.997 (0.997 - 0.998)	0.360 (0.334 - 0.386)	0.429 (0.405 - 0.451)
	Merl-R18	0.659 (0.635 - 0.682)	0.926 (0.919 - 0.933)	0.220 (0.199 - 0.240)	0.998 (0.998 - 0.999)	0.354 (0.327 - 0.380)	0.430 (0.406 - 0.452)
	Random Init-R18	0.921 (0.903 - 0.938)	0.986 (0.979 - 0.992)	0.762 (0.729 - 0.794)	0.999 (0.998 - 0.999)	0.847 (0.825 - 0.869)	0.848 (0.827 - 0.870)
	ECGCLIP-R18	0.923 (0.905 - 0.939)	0.985 (0.977 - 0.991)	0.752 (0.719 - 0.788)	0.999 (0.999 - 1.000)	0.847 (0.823 - 0.870)	0.850 (0.828 - 0.871)
	ECGCLIP-R34	0.925 (0.908 - 0.940)	0.984 (0.977 - 0.991)	0.778 (0.745 - 0.810)	0.999 (0.999 - 0.999)	0.859 (0.837 - 0.879)	0.860 (0.839 - 0.880)
	Merl-R18	0.920 (0.901 - 0.936)	0.985 (0.977 - 0.991)	0.732 (0.697 - 0.766)	0.999 (0.999 - 1.000)	0.833 (0.809 - 0.856)	0.837 (0.814 - 0.858)
IVB	Random Init-R18	0.122 (0.107 - 0.141)	0.728 (0.707 - 0.749)	0.080 (0.062 - 0.099)	0.989 (0.987 - 0.990)	0.116 (0.091 - 0.140)	0.110 (0.083 - 0.138)
	ECGCLIP-R18	0.133 (0.115 - 0.153)	0.742 (0.722 - 0.761)	0.085 (0.066 - 0.104)	0.989 (0.988 - 0.991)	0.125 (0.097 - 0.151)	0.122 (0.091 - 0.151)
	ECGCLIP-R34	0.125 (0.108 - 0.145)	0.726 (0.705 - 0.745)	0.079 (0.059 - 0.097)	0.991 (0.989 - 0.992)	0.119 (0.090 - 0.145)	0.119 (0.087 - 0.149)
	Merl-R18	0.121 (0.105 - 0.138)	0.732 (0.710 - 0.751)	0.079 (0.060 - 0.096)	0.989 (0.988 - 0.991)	0.115 (0.089 - 0.140)	0.111 (0.082 - 0.138)
VPE	Random Init-R18	0.729 (0.621 - 0.828)	0.962 (0.927 - 0.990)	0.544 (0.432 - 0.659)	1.000 (1.000 - 1.000)	0.688 (0.586 - 0.776)	0.713 (0.628 - 0.789)
	ECGCLIP-R18	0.758 (0.662 - 0.849)	0.960 (0.928 - 0.987)	0.633 (0.521 - 0.738)	1.000 (0.999 - 1.000)	0.741 (0.649 - 0.819)	0.751 (0.666 - 0.826)
	ECGCLIP-R34	0.761 (0.664 - 0.853)	0.963 (0.929 - 0.990)	0.608 (0.493 - 0.710)	1.000 (1.000 - 1.000)	0.744 (0.646 - 0.819)	0.763 (0.684 - 0.832)
	Merl-R18	0.704 (0.596 - 0.809)	0.961 (0.927 - 0.988)	0.646 (0.530 - 0.754)	1.000 (0.999 - 1.000)	0.729 (0.639 - 0.805)	0.734 (0.648 - 0.809)
3° AVB	Random Init-R18	0.135 (0.038 - 0.314)	0.989 (0.981 - 0.996)	0.500 (0.250 - 0.750)	0.997 (0.996 - 0.998)	0.168 (0.069 - 0.278)	0.224 (0.104 - 0.345)
	ECGCLIP-R18	0.361 (0.101 - 0.598)	0.993 (0.988 - 0.998)	0.625 (0.375 - 0.875)	0.996 (0.995 - 0.996)	0.161 (0.073 - 0.248)	0.239 (0.128 - 0.339)
	ECGCLIP-R34	0.502 (0.235 - 0.732)	0.994 (0.989 - 0.998)	0.500 (0.250 - 0.739)	1.000 (1.000 - 1.000)	0.552 (0.308 - 0.759)	0.554 (0.311 - 0.763)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
2° 1 Type AVB	Merl-R18	0.053 (0.017 - 0.203)	0.987 (0.980 - 0.993)	0.375 (0.143 - 0.611)	0.996 (0.995 - 0.997)	0.106 (0.035 - 0.188)	0.151 (0.053 - 0.250)
	Random Init-R18	0.169 (0.010 - 0.402)	0.788 (0.644 - 0.922)	0.214 (0.000 - 0.429)	0.999 (0.999 - 1.000)	0.214 (0.000 - 0.390)	0.214 (0.000 - 0.398)
	ECGCLIP-R18	0.076 (0.008 - 0.255)	0.835 (0.697 - 0.949)	0.286 (0.063 - 0.500)	0.998 (0.997 - 0.999)	0.131 (0.031 - 0.238)	0.155 (0.036 - 0.272)
	ECGCLIP-R34	0.053 (0.006 - 0.176)	0.824 (0.687 - 0.942)	0.286 (0.063 - 0.500)	0.999 (0.998 - 0.999)	0.160 (0.038 - 0.292)	0.177 (0.043 - 0.314)
2° 2 Type AVB	Merl-R18	0.057 (0.006 - 0.210)	0.794 (0.630 - 0.936)	0.214 (0.000 - 0.429)	0.999 (0.999 - 1.000)	0.194 (0.000 - 0.357)	0.194 (-0.001 - 0.365)
	Random Init-R18	0.151 (0.002 - 0.368)	0.726 (0.551 - 0.882)	0.143 (0.000 - 0.357)	1.000 (1.000 - 1.000)	0.200 (0.000 - 0.435)	0.218 (0.000 - 0.456)
	ECGCLIP-R18	0.058 (0.005 - 0.268)	0.846 (0.721 - 0.943)	0.143 (0.000 - 0.357)	1.000 (1.000 - 1.000)	0.190 (0.000 - 0.417)	0.202 (0.000 - 0.436)
	ECGCLIP-R34	0.042 (0.005 - 0.193)	0.634 (0.413 - 0.839)	0.071 (0.000 - 0.250)	1.000 (0.999 - 1.000)	0.091 (0.000 - 0.273)	0.094 (-0.001 - 0.283)
LVH	Merl-R18	0.155 (0.003 - 0.373)	0.656 (0.446 - 0.851)	0.143 (0.000 - 0.357)	1.000 (1.000 - 1.000)	0.222 (0.000 - 0.476)	0.267 (0.000 - 0.527)
	Random Init-R18	0.649 (0.630 - 0.667)	0.891 (0.883 - 0.898)	0.378 (0.358 - 0.397)	0.988 (0.987 - 0.990)	0.513 (0.493 - 0.533)	0.516 (0.497 - 0.534)
	ECGCLIP-R18	0.660 (0.642 - 0.680)	0.899 (0.892 - 0.906)	0.382 (0.363 - 0.401)	0.989 (0.988 - 0.991)	0.519 (0.500 - 0.539)	0.524 (0.505 - 0.544)
	ECGCLIP-R34	0.643 (0.623 - 0.663)	0.893 (0.886 - 0.901)	0.433 (0.413 - 0.452)	0.983 (0.981 - 0.985)	0.550 (0.530 - 0.568)	0.535 (0.515 - 0.553)
LAH	Merl-R18	0.639 (0.620 - 0.658)	0.886 (0.878 - 0.894)	0.411 (0.390 - 0.430)	0.985 (0.983 - 0.987)	0.536 (0.514 - 0.555)	0.527 (0.507 - 0.546)
	Random Init-R18	0.147 (0.121 - 0.175)	0.871 (0.855 - 0.885)	0.333 (0.288 - 0.379)	0.964 (0.961 - 0.966)	0.211 (0.183 - 0.240)	0.205 (0.175 - 0.236)
	ECGCLIP-R18	0.140 (0.117 - 0.169)	0.867 (0.851 - 0.882)	0.317 (0.274 - 0.363)	0.967 (0.965 - 0.970)	0.215 (0.183 - 0.244)	0.205 (0.173 - 0.237)
	ECGCLIP-R34	0.146 (0.122 - 0.174)	0.874 (0.861 - 0.888)	0.357 (0.311 - 0.402)	0.959 (0.956 - 0.961)	0.208 (0.180 - 0.236)	0.205 (0.174 - 0.235)
	Merl-R18	0.138 (0.115 - 0.166)	0.864 (0.848 - 0.879)	0.378 (0.333 - 0.424)	0.956 (0.953 - 0.959)	0.211 (0.186 - 0.238)	0.211 (0.184 - 0.240)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RVH	Random Init-R18	0.107 (0.061 – 0.168)	0.918 (0.899 – 0.935)	0.151 (0.092 – 0.219)	0.996 (0.995 – 0.997)	0.160 (0.098 – 0.229)	0.155 (0.093 – 0.225)
	ECGCLIP-R18	0.112 (0.067 – 0.180)	0.907 (0.885 – 0.927)	0.127 (0.071 – 0.189)	0.997 (0.996 – 0.998)	0.158 (0.092 – 0.229)	0.159 (0.089 – 0.231)
	ECGCLIP-R34	0.090 (0.052 – 0.147)	0.906 (0.885 – 0.927)	0.111 (0.058 – 0.169)	0.997 (0.996 – 0.998)	0.138 (0.075 – 0.205)	0.138 (0.075 – 0.207)
	Merl-R18	0.109 (0.064 – 0.174)	0.924 (0.905 – 0.940)	0.119 (0.066 – 0.180)	0.998 (0.997 – 0.998)	0.160 (0.091 – 0.236)	0.166 (0.095 – 0.247)
VA	Random Init-R18	0.054 (0.037 – 0.077)	0.910 (0.882 – 0.934)	0.144 (0.074 – 0.217)	0.992 (0.991 – 0.993)	0.106 (0.056 – 0.158)	0.104 (0.051 – 0.160)
	ECGCLIP-R18	0.045 (0.032 – 0.064)	0.914 (0.892 – 0.934)	0.067 (0.022 – 0.118)	0.994 (0.993 – 0.995)	0.058 (0.018 – 0.099)	0.053 (0.013 – 0.096)
	ECGCLIP-R34	0.060 (0.041 – 0.086)	0.924 (0.904 – 0.941)	0.163 (0.084 – 0.239)	0.994 (0.993 – 0.995)	0.133 (0.071 – 0.190)	0.131 (0.067 – 0.188)
	Merl-R18	0.047 (0.033 – 0.070)	0.920 (0.901 – 0.937)	0.144 (0.080 – 0.215)	0.991 (0.990 – 0.992)	0.095 (0.053 – 0.142)	0.095 (0.050 – 0.145)
RAH	Random Init-R18	0.109 (0.057 – 0.175)	0.903 (0.868 – 0.931)	0.101 (0.044 – 0.161)	1.000 (0.999 – 1.000)	0.171 (0.077 – 0.258)	0.236 (0.124 – 0.337)
	ECGCLIP-R18	0.163 (0.097 – 0.244)	0.929 (0.895 – 0.954)	0.101 (0.045 – 0.167)	1.000 (0.999 – 1.000)	0.169 (0.079 – 0.264)	0.229 (0.119 – 0.328)
	ECGCLIP-R34	0.217 (0.140 – 0.304)	0.925 (0.889 – 0.954)	0.051 (0.011 – 0.095)	1.000 (1.000 – 1.000)	0.093 (0.021 – 0.170)	0.167 (0.048 – 0.270)
	Merl-R18	0.144 (0.081 – 0.225)	0.924 (0.900 – 0.944)	0.081 (0.033 – 0.134)	1.000 (1.000 – 1.000)	0.145 (0.063 – 0.232)	0.242 (0.131 – 0.343)
VVI	Random Init-R18	0.761 (0.711 – 0.807)	0.966 (0.952 – 0.978)	0.667 (0.609 – 0.718)	0.998 (0.997 – 0.998)	0.730 (0.687 – 0.771)	0.730 (0.687 – 0.773)
	ECGCLIP-R18	0.763 (0.716 – 0.808)	0.968 (0.955 – 0.980)	0.660 (0.607 – 0.714)	0.998 (0.998 – 0.999)	0.736 (0.692 – 0.780)	0.738 (0.695 – 0.782)
	ECGCLIP-R34	0.701 (0.653 – 0.745)	0.948 (0.931 – 0.963)	0.595 (0.540 – 0.646)	0.998 (0.997 – 0.999)	0.686 (0.637 – 0.730)	0.691 (0.644 – 0.734)
	Merl-R18	0.784 (0.737 – 0.827)	0.971 (0.958 – 0.982)	0.626 (0.566 – 0.679)	0.999 (0.999 – 0.999)	0.737 (0.692 – 0.778)	0.747 (0.705 – 0.785)
VAT	Random Init-R18	0.472 (0.410 – 0.532)	0.927 (0.905 – 0.946)	0.160 (0.119 – 0.204)	0.999 (0.999 – 1.000)	0.266 (0.205 – 0.326)	0.351 (0.285 – 0.409)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
DDD	ECGCLIP-R18	0.504 (0.445 – 0.565)	0.935 (0.916 – 0.952)	0.129 (0.092 – 0.171)	1.000 (0.999 – 1.000)	0.223 (0.164 – 0.287)	0.320 (0.255 – 0.382)
	ECGCLIP-R34	0.426 (0.366 – 0.488)	0.909 (0.883 – 0.932)	0.109 (0.074 – 0.143)	1.000 (0.999 – 1.000)	0.191 (0.134 – 0.245)	0.289 (0.223 – 0.347)
	Merl-R18	0.557 (0.493 – 0.618)	0.932 (0.910 – 0.951)	0.180 (0.138 – 0.225)	0.999 (0.999 – 1.000)	0.295 (0.233 – 0.356)	0.380 (0.320 – 0.437)
	Random Init-R18	0.761 (0.711 – 0.809)	0.958 (0.941 – 0.973)	0.163 (0.123 – 0.211)	1.000 (1.000 – 1.000)	0.278 (0.218 – 0.346)	0.389 (0.335 – 0.448)
	ECGCLIP-R18	0.652 (0.592 – 0.709)	0.960 (0.944 – 0.975)	0.160 (0.118 – 0.202)	1.000 (1.000 – 1.000)	0.272 (0.208 – 0.331)	0.377 (0.321 – 0.431)
	ECGCLIP-R34	0.588 (0.526 – 0.647)	0.965 (0.952 – 0.977)	0.126 (0.090 – 0.164)	0.999 (0.999 – 1.000)	0.216 (0.160 – 0.270)	0.308 (0.246 – 0.364)
	Merl-R18	0.792 (0.743 – 0.839)	0.966 (0.950 – 0.980)	0.088 (0.057 – 0.122)	1.000 (1.000 – 1.000)	0.162 (0.108 – 0.217)	0.290 (0.231 – 0.344)
	Random Init-R18	0.152 (0.110 – 0.199)	0.813 (0.786 – 0.840)	0.088 (0.057 – 0.125)	0.999 (0.999 – 0.999)	0.153 (0.101 – 0.209)	0.220 (0.156 – 0.283)
	ECGCLIP-R18	0.053 (0.038 – 0.077)	0.720 (0.691 – 0.750)	0.048 (0.025 – 0.075)	0.998 (0.997 – 0.998)	0.078 (0.042 – 0.121)	0.097 (0.051 – 0.149)
	ECGCLIP-R34	0.083 (0.051 – 0.127)	0.345 (0.306 – 0.386)	0.085 (0.054 – 0.122)	0.999 (0.998 – 0.999)	0.144 (0.093 – 0.200)	0.194 (0.134 – 0.259)
	Merl-R18	0.071 (0.051 – 0.103)	0.747 (0.718 – 0.775)	0.024 (0.007 – 0.045)	0.999 (0.999 – 1.000)	0.045 (0.014 – 0.082)	0.094 (0.036 – 0.161)
	Random Init-R18	0.782 (0.739 – 0.824)	0.969 (0.956 – 0.981)	0.636 (0.580 – 0.686)	0.999 (0.999 – 1.000)	0.751 (0.706 – 0.790)	0.761 (0.719 – 0.798)
VP	ECGCLIP-R18	0.770 (0.728 – 0.810)	0.971 (0.958 – 0.982)	0.558 (0.504 – 0.611)	1.000 (0.999 – 1.000)	0.701 (0.653 – 0.744)	0.722 (0.681 – 0.759)
	ECGCLIP-R34	0.737 (0.692 – 0.779)	0.962 (0.947 – 0.976)	0.265 (0.217 – 0.318)	1.000 (1.000 – 1.000)	0.419 (0.357 – 0.483)	0.513 (0.464 – 0.561)
	Merl-R18	0.805 (0.765 – 0.844)	0.972 (0.958 – 0.983)	0.636 (0.582 – 0.689)	0.999 (0.999 – 1.000)	0.751 (0.708 – 0.790)	0.761 (0.722 – 0.797)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S24: Detailed task-specific performance comparison on the external Georgia cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STEMI	Random Init-R18	0.002 (0.000 - 0.004)	0.707 (0.396 - 0.917)	0.000 (0.000 - 0.000)	0.996 (0.994 - 0.997)	0.000 (0.000 - 0.000)	-0.002 (-0.002 - -0.001)
	ECGCLIP-R18	0.002 (0.000 - 0.006)	0.717 (0.411 - 0.919)	0.000 (0.000 - 0.000)	0.995 (0.993 - 0.996)	0.000 (0.000 - 0.000)	-0.002 (-0.003 - -0.001)
	ECGCLIP-R34	0.004 (0.000 - 0.025)	0.698 (0.382 - 0.909)	0.143 (0.000 - 0.500)	0.991 (0.990 - 0.993)	0.020 (0.000 - 0.065)	0.037 (-0.003 - 0.119)
	Merl-R18	0.002 (0.000 - 0.005)	0.736 (0.490 - 0.900)	0.000 (0.000 - 0.000)	0.994 (0.992 - 0.995)	0.000 (0.000 - 0.000)	-0.002 (-0.003 - -0.001)
SBrad	Random Init-R18	0.934 (0.923 - 0.944)	0.983 (0.980 - 0.987)	0.872 (0.855 - 0.888)	0.975 (0.972 - 0.979)	0.872 (0.860 - 0.884)	0.848 (0.833 - 0.862)
	ECGCLIP-R18	0.939 (0.928 - 0.949)	0.985 (0.982 - 0.988)	0.859 (0.842 - 0.875)	0.979 (0.976 - 0.982)	0.874 (0.861 - 0.886)	0.850 (0.835 - 0.864)
	ECGCLIP-R34	0.967 (0.959 - 0.974)	0.992 (0.990 - 0.994)	0.924 (0.912 - 0.936)	0.985 (0.983 - 0.988)	0.925 (0.915 - 0.934)	0.910 (0.899 - 0.921)
	Merl-R18	0.928 (0.917 - 0.939)	0.981 (0.978 - 0.985)	0.834 (0.815 - 0.851)	0.981 (0.978 - 0.984)	0.863 (0.850 - 0.875)	0.838 (0.824 - 0.853)
STach	Random Init-R18	0.934 (0.924 - 0.944)	0.988 (0.985 - 0.990)	0.831 (0.812 - 0.851)	0.988 (0.986 - 0.990)	0.867 (0.854 - 0.881)	0.850 (0.836 - 0.866)
	ECGCLIP-R18	0.936 (0.926 - 0.947)	0.987 (0.984 - 0.990)	0.830 (0.810 - 0.851)	0.989 (0.987 - 0.991)	0.868 (0.855 - 0.883)	0.852 (0.837 - 0.868)
	ECGCLIP-R34	0.956 (0.946 - 0.965)	0.993 (0.991 - 0.994)	0.880 (0.864 - 0.897)	0.992 (0.991 - 0.994)	0.909 (0.898 - 0.921)	0.898 (0.886 - 0.911)
	Merl-R18	0.932 (0.921 - 0.943)	0.987 (0.983 - 0.989)	0.838 (0.819 - 0.859)	0.987 (0.985 - 0.989)	0.867 (0.854 - 0.882)	0.850 (0.834 - 0.867)
VPB	Random Init-R18	0.358 (0.305 - 0.406)	0.810 (0.785 - 0.832)	0.371 (0.321 - 0.418)	0.974 (0.971 - 0.977)	0.367 (0.320 - 0.407)	0.342 (0.295 - 0.383)
	ECGCLIP-R18	0.370 (0.319 - 0.423)	0.814 (0.790 - 0.836)	0.371 (0.324 - 0.418)	0.978 (0.975 - 0.981)	0.386 (0.340 - 0.429)	0.363 (0.317 - 0.408)
	ECGCLIP-R34	0.441 (0.388 - 0.489)	0.835 (0.811 - 0.858)	0.421 (0.371 - 0.471)	0.983 (0.980 - 0.985)	0.455 (0.410 - 0.499)	0.437 (0.391 - 0.482)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
APB	Merl-R18	0.352 (0.302 - 0.402)	0.808 (0.784 - 0.831)	0.373 (0.325 - 0.421)	0.979 (0.977 - 0.982)	0.396 (0.349 - 0.439)	0.374 (0.327 - 0.419)
	Random Init-R18	0.146 (0.129 - 0.164)	0.750 (0.733 - 0.767)	0.538 (0.502 - 0.578)	0.784 (0.775 - 0.792)	0.223 (0.204 - 0.243)	0.182 (0.161 - 0.206)
	ECGCLIP-R18	0.180 (0.158 - 0.205)	0.779 (0.761 - 0.795)	0.481 (0.444 - 0.522)	0.852 (0.845 - 0.860)	0.258 (0.236 - 0.281)	0.214 (0.189 - 0.240)
	ECGCLIP-R34	0.305 (0.271 - 0.342)	0.837 (0.822 - 0.854)	0.487 (0.448 - 0.530)	0.923 (0.917 - 0.928)	0.366 (0.336 - 0.395)	0.325 (0.294 - 0.355)
SArr	Merl-R18	0.157 (0.138 - 0.177)	0.762 (0.745 - 0.779)	0.520 (0.481 - 0.560)	0.816 (0.808 - 0.824)	0.241 (0.219 - 0.260)	0.201 (0.175 - 0.224)
	Random Init-R18	0.096 (0.082 - 0.116)	0.719 (0.699 - 0.740)	0.037 (0.021 - 0.054)	0.992 (0.991 - 0.994)	0.062 (0.036 - 0.090)	0.065 (0.029 - 0.103)
	ECGCLIP-R18	0.100 (0.085 - 0.123)	0.723 (0.703 - 0.744)	0.066 (0.043 - 0.090)	0.987 (0.985 - 0.989)	0.098 (0.064 - 0.130)	0.090 (0.052 - 0.128)
	ECGCLIP-R34	0.110 (0.092 - 0.134)	0.735 (0.716 - 0.755)	0.104 (0.075 - 0.132)	0.971 (0.968 - 0.974)	0.119 (0.088 - 0.150)	0.086 (0.054 - 0.119)
AF	Merl-R18	0.100 (0.082 - 0.121)	0.707 (0.686 - 0.729)	0.073 (0.049 - 0.097)	0.985 (0.983 - 0.988)	0.104 (0.071 - 0.137)	0.091 (0.055 - 0.128)
	Random Init-R18	0.694 (0.646 - 0.737)	0.922 (0.906 - 0.937)	0.595 (0.557 - 0.632)	0.993 (0.991 - 0.994)	0.691 (0.659 - 0.721)	0.686 (0.653 - 0.716)
	ECGCLIP-R18	0.710 (0.665 - 0.752)	0.922 (0.905 - 0.937)	0.574 (0.531 - 0.614)	0.995 (0.993 - 0.996)	0.691 (0.655 - 0.723)	0.692 (0.660 - 0.725)
	ECGCLIP-R34	0.746 (0.707 - 0.782)	0.927 (0.911 - 0.941)	0.626 (0.588 - 0.663)	0.995 (0.993 - 0.996)	0.729 (0.697 - 0.758)	0.727 (0.696 - 0.755)
AFL	Merl-R18	0.701 (0.655 - 0.744)	0.921 (0.905 - 0.936)	0.611 (0.571 - 0.649)	0.993 (0.992 - 0.995)	0.709 (0.677 - 0.740)	0.705 (0.674 - 0.736)
	Random Init-R18	0.619 (0.537 - 0.692)	0.935 (0.907 - 0.958)	0.575 (0.500 - 0.649)	0.995 (0.994 - 0.996)	0.622 (0.557 - 0.681)	0.618 (0.553 - 0.678)
	ECGCLIP-R18	0.672 (0.604 - 0.736)	0.938 (0.911 - 0.961)	0.640 (0.574 - 0.707)	0.994 (0.992 - 0.995)	0.645 (0.587 - 0.702)	0.639 (0.581 - 0.696)
	ECGCLIP-R34	0.656 (0.584 - 0.728)	0.943 (0.916 - 0.965)	0.618 (0.549 - 0.685)	0.994 (0.992 - 0.995)	0.630 (0.571 - 0.688)	0.624 (0.564 - 0.683)
	Merl-R18	0.637 (0.567 - 0.708)	0.935 (0.907 - 0.959)	0.586 (0.518 - 0.656)	0.995 (0.993 - 0.996)	0.626 (0.567 - 0.683)	0.622 (0.562 - 0.679)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
JPB	Random Init-R18	0.000 (0.000 - 0.001)	0.750 (0.742 - 0.759)	0.000 (0.000 - 0.000)	0.909 (0.904 - 0.915)	0.000 (0.000 - 0.000)	-0.003 (-0.006 - -0.003)
	ECGCLIP-R18	0.000 (0.000 - 0.001)	0.821 (0.815 - 0.829)	0.000 (0.000 - 0.000)	0.957 (0.954 - 0.961)	0.000 (0.000 - 0.000)	-0.002 (-0.004 - -0.002)
	ECGCLIP-R34	0.001 (0.001 - 0.002)	0.925 (0.920 - 0.930)	0.000 (0.000 - 0.000)	0.981 (0.978 - 0.983)	0.000 (0.000 - 0.000)	-0.001 (-0.003 - -0.001)
	Merl-R18	0.000 (0.000 - 0.001)	0.774 (0.766 - 0.782)	0.000 (0.000 - 0.000)	0.975 (0.971 - 0.977)	0.000 (0.000 - 0.000)	-0.002 (-0.003 - -0.002)
AJR	Random Init-R18	0.009 (0.004 - 0.016)	0.856 (0.764 - 0.923)	0.053 (0.000 - 0.182)	0.981 (0.979 - 0.984)	0.009 (0.000 - 0.031)	0.011 (-0.007 - 0.050)
	ECGCLIP-R18	0.011 (0.005 - 0.024)	0.860 (0.779 - 0.924)	0.211 (0.048 - 0.412)	0.979 (0.976 - 0.982)	0.033 (0.008 - 0.066)	0.056 (0.008 - 0.113)
	ECGCLIP-R34	0.009 (0.004 - 0.018)	0.855 (0.785 - 0.916)	0.053 (0.000 - 0.182)	0.990 (0.988 - 0.992)	0.016 (0.000 - 0.052)	0.018 (-0.005 - 0.067)
	Merl-R18	0.013 (0.006 - 0.027)	0.882 (0.803 - 0.937)	0.158 (0.000 - 0.333)	0.987 (0.984 - 0.989)	0.038 (0.000 - 0.082)	0.054 (-0.005 - 0.119)
ATach	Random Init-R18	0.031 (0.014 - 0.068)	0.887 (0.807 - 0.952)	0.571 (0.364 - 0.750)	0.958 (0.954 - 0.961)	0.066 (0.034 - 0.099)	0.134 (0.076 - 0.184)
	ECGCLIP-R18	0.041 (0.017 - 0.089)	0.894 (0.826 - 0.957)	0.429 (0.250 - 0.621)	0.975 (0.972 - 0.977)	0.079 (0.040 - 0.122)	0.130 (0.066 - 0.194)
	ECGCLIP-R34	0.051 (0.025 - 0.091)	0.927 (0.863 - 0.975)	0.571 (0.381 - 0.750)	0.974 (0.971 - 0.977)	0.102 (0.057 - 0.147)	0.172 (0.106 - 0.236)
	Merl-R18	0.033 (0.015 - 0.062)	0.913 (0.864 - 0.957)	0.429 (0.240 - 0.619)	0.972 (0.969 - 0.975)	0.072 (0.033 - 0.113)	0.123 (0.060 - 0.182)
JTach	Random Init-R18	0.006 (0.000 - 0.033)	0.790 (0.315 - 0.995)	0.250 (0.000 - 1.000)	0.994 (0.993 - 0.996)	0.030 (0.000 - 0.107)	0.062 (-0.002 - 0.188)
	ECGCLIP-R18	0.003 (0.000 - 0.016)	0.813 (0.456 - 0.990)	0.000 (0.000 - 0.000)	0.994 (0.992 - 0.996)	0.000 (0.000 - 0.000)	-0.002 (-0.002 - -0.001)
	ECGCLIP-R34	0.011 (0.000 - 0.070)	0.810 (0.536 - 0.998)	0.250 (0.000 - 1.000)	0.995 (0.993 - 0.996)	0.033 (0.000 - 0.118)	0.065 (-0.002 - 0.201)
	Merl-R18	0.007 (0.000 - 0.038)	0.782 (0.217 - 0.996)	0.250 (0.000 - 1.000)	0.993 (0.991 - 0.995)	0.026 (0.000 - 0.092)	0.057 (-0.002 - 0.173)
JEB	Random Init-R18	0.019 (0.000 - 0.083)	0.765 (0.439 - 0.998)	0.400 (0.000 - 1.000)	0.992 (0.990 - 0.994)	0.044 (0.000 - 0.109)	0.095 (-0.002 - 0.212)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
SVT	ECGCLIP-R18	0.014 (0.000 - 0.062)	0.824 (0.621 - 0.997)	0.400 (0.000 - 1.000)	0.994 (0.993 - 0.996)	0.063 (0.000 - 0.152)	0.115 (-0.002 - 0.256)
	ECGCLIP-R34	0.034 (0.000 - 0.158)	0.849 (0.642 - 0.999)	0.400 (0.000 - 1.000)	0.997 (0.996 - 0.998)	0.103 (0.000 - 0.244)	0.152 (-0.001 - 0.329)
	Merl-R18	0.018 (0.000 - 0.079)	0.802 (0.564 - 0.998)	0.400 (0.000 - 1.000)	0.997 (0.995 - 0.998)	0.093 (0.000 - 0.227)	0.144 (-0.001 - 0.310)
	Random Init-R18	0.324 (0.174 - 0.505)	0.947 (0.894 - 0.987)	0.250 (0.103 - 0.412)	0.999 (0.999 - 1.000)	0.327 (0.154 - 0.491)	0.342 (0.160 - 0.505)
	ECGCLIP-R18	0.348 (0.174 - 0.523)	0.943 (0.890 - 0.987)	0.156 (0.037 - 0.294)	1.000 (0.999 - 1.000)	0.238 (0.057 - 0.409)	0.278 (0.077 - 0.459)
	ECGCLIP-R34	0.331 (0.175 - 0.506)	0.940 (0.878 - 0.989)	0.188 (0.064 - 0.333)	0.999 (0.999 - 1.000)	0.261 (0.093 - 0.429)	0.282 (0.106 - 0.455)
	Merl-R18	0.313 (0.160 - 0.507)	0.919 (0.851 - 0.973)	0.188 (0.057 - 0.324)	0.999 (0.999 - 1.000)	0.267 (0.089 - 0.429)	0.293 (0.107 - 0.464)
	Random Init-R18	0.000 (0.000 - 0.001)	0.772 (0.764 - 0.780)	0.000 (0.000 - 0.000)	0.999 (0.998 - 0.999)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
VEB	ECGCLIP-R18	0.000 (0.000 - 0.000)	0.522 (0.513 - 0.532)	0.000 (0.000 - 0.000)	0.999 (0.998 - 1.000)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
	ECGCLIP-R34	0.000 (0.000 - 0.001)	0.647 (0.637 - 0.656)	0.000 (0.000 - 0.000)	0.999 (0.998 - 0.999)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
	Merl-R18	0.000 (0.000 - 0.001)	0.757 (0.749 - 0.765)	0.000 (0.000 - 0.000)	0.995 (0.994 - 0.996)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)
	Random Init-R18	0.794 (0.768 - 0.820)	0.939 (0.930 - 0.949)	0.748 (0.722 - 0.775)	0.974 (0.971 - 0.977)	0.751 (0.731 - 0.772)	0.725 (0.704 - 0.747)
RBBB	ECGCLIP-R18	0.792 (0.768 - 0.817)	0.942 (0.934 - 0.951)	0.722 (0.697 - 0.751)	0.976 (0.973 - 0.980)	0.742 (0.721 - 0.764)	0.716 (0.693 - 0.739)
	ECGCLIP-R34	0.794 (0.770 - 0.819)	0.947 (0.938 - 0.956)	0.729 (0.703 - 0.758)	0.973 (0.969 - 0.976)	0.733 (0.713 - 0.754)	0.705 (0.684 - 0.729)
	Merl-R18	0.809 (0.785 - 0.832)	0.953 (0.946 - 0.961)	0.743 (0.718 - 0.769)	0.975 (0.972 - 0.979)	0.753 (0.732 - 0.773)	0.727 (0.705 - 0.749)
1° AVB	Random Init-R18	0.758 (0.725 - 0.791)	0.935 (0.922 - 0.947)	0.705 (0.674 - 0.733)	0.981 (0.978 - 0.983)	0.732 (0.708 - 0.754)	0.710 (0.684 - 0.734)
	ECGCLIP-R18	0.758 (0.725 - 0.793)	0.926 (0.912 - 0.940)	0.726 (0.694 - 0.755)	0.979 (0.976 - 0.982)	0.741 (0.715 - 0.763)	0.719 (0.691 - 0.742)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
QT Prolong	ECGCLIP-R34	0.767 (0.735 - 0.800)	0.927 (0.913 - 0.941)	0.698 (0.666 - 0.726)	0.985 (0.982 - 0.987)	0.747 (0.722 - 0.769)	0.729 (0.703 - 0.752)
	Merl-R18	0.766 (0.734 - 0.797)	0.934 (0.922 - 0.946)	0.746 (0.718 - 0.773)	0.977 (0.974 - 0.980)	0.744 (0.719 - 0.766)	0.721 (0.695 - 0.745)
	Random Init-R18	0.508 (0.480 - 0.535)	0.844 (0.833 - 0.854)	0.224 (0.203 - 0.248)	0.982 (0.979 - 0.985)	0.334 (0.309 - 0.364)	0.337 (0.311 - 0.368)
	ECGCLIP-R18	0.551 (0.525 - 0.578)	0.863 (0.853 - 0.872)	0.261 (0.239 - 0.286)	0.984 (0.981 - 0.986)	0.383 (0.356 - 0.412)	0.388 (0.360 - 0.417)
	ECGCLIP-R34	0.613 (0.587 - 0.638)	0.887 (0.878 - 0.896)	0.386 (0.360 - 0.410)	0.977 (0.973 - 0.980)	0.502 (0.477 - 0.526)	0.477 (0.452 - 0.502)
	Merl-R18	0.510 (0.483 - 0.538)	0.843 (0.832 - 0.853)	0.222 (0.201 - 0.244)	0.985 (0.982 - 0.987)	0.336 (0.309 - 0.365)	0.347 (0.318 - 0.377)
ER	Random Init-R18	0.305 (0.218 - 0.396)	0.905 (0.875 - 0.932)	0.129 (0.076 - 0.188)	0.998 (0.998 - 0.999)	0.207 (0.125 - 0.292)	0.256 (0.160 - 0.350)
	ECGCLIP-R18	0.328 (0.243 - 0.414)	0.924 (0.898 - 0.946)	0.057 (0.021 - 0.094)	0.999 (0.999 - 1.000)	0.103 (0.039 - 0.163)	0.166 (0.070 - 0.248)
	ECGCLIP-R34	0.316 (0.240 - 0.401)	0.927 (0.900 - 0.949)	0.186 (0.124 - 0.248)	0.998 (0.997 - 0.999)	0.275 (0.190 - 0.351)	0.309 (0.220 - 0.389)
	Merl-R18	0.281 (0.195 - 0.363)	0.880 (0.842 - 0.914)	0.121 (0.069 - 0.178)	0.999 (0.998 - 0.999)	0.199 (0.120 - 0.278)	0.254 (0.163 - 0.338)
LAFB	Random Init-R18	0.205 (0.164 - 0.256)	0.934 (0.917 - 0.947)	0.289 (0.225 - 0.359)	0.987 (0.985 - 0.989)	0.283 (0.224 - 0.343)	0.270 (0.211 - 0.331)
	ECGCLIP-R18	0.221 (0.179 - 0.276)	0.939 (0.923 - 0.953)	0.194 (0.140 - 0.252)	0.993 (0.991 - 0.994)	0.243 (0.178 - 0.304)	0.241 (0.173 - 0.303)
	ECGCLIP-R34	0.268 (0.214 - 0.334)	0.945 (0.929 - 0.958)	0.250 (0.189 - 0.314)	0.991 (0.990 - 0.993)	0.288 (0.224 - 0.351)	0.280 (0.218 - 0.347)
LBBB	Merl-R18	0.212 (0.168 - 0.268)	0.931 (0.914 - 0.945)	0.206 (0.147 - 0.267)	0.991 (0.989 - 0.993)	0.240 (0.176 - 0.303)	0.233 (0.169 - 0.296)
	Random Init-R18	0.674 (0.622 - 0.726)	0.943 (0.927 - 0.959)	0.373 (0.323 - 0.427)	0.998 (0.997 - 0.999)	0.524 (0.472 - 0.578)	0.564 (0.519 - 0.612)
	ECGCLIP-R18	0.690 (0.638 - 0.740)	0.940 (0.923 - 0.955)	0.367 (0.319 - 0.420)	0.999 (0.998 - 1.000)	0.525 (0.472 - 0.582)	0.575 (0.531 - 0.623)
	ECGCLIP-R34	0.661 (0.605 - 0.714)	0.932 (0.912 - 0.950)	0.385 (0.335 - 0.440)	0.998 (0.997 - 0.999)	0.536 (0.482 - 0.591)	0.575 (0.527 - 0.623)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
IVB	Merl-R18	0.678 (0.621 - 0.733)	0.932 (0.911 - 0.950)	0.361 (0.315 - 0.417)	0.999 (0.998 - 1.000)	0.518 (0.467 - 0.574)	0.567 (0.523 - 0.614)
	Random Init-R18	0.066 (0.055 - 0.085)	0.807 (0.780 - 0.835)	0.084 (0.048 - 0.126)	0.984 (0.981 - 0.986)	0.088 (0.051 - 0.130)	0.071 (0.033 - 0.113)
	ECGCLIP-R18	0.066 (0.053 - 0.084)	0.801 (0.773 - 0.830)	0.064 (0.032 - 0.100)	0.983 (0.981 - 0.985)	0.067 (0.034 - 0.103)	0.049 (0.017 - 0.086)
	ECGCLIP-R34	0.056 (0.046 - 0.071)	0.766 (0.733 - 0.798)	0.039 (0.015 - 0.069)	0.984 (0.981 - 0.986)	0.043 (0.016 - 0.073)	0.025 (-0.001 - 0.055)
Short PR	Merl-R18	0.063 (0.052 - 0.081)	0.801 (0.773 - 0.829)	0.044 (0.019 - 0.073)	0.984 (0.982 - 0.986)	0.048 (0.021 - 0.078)	0.031 (0.003 - 0.061)
	Random Init-R18	0.002 (0.001 - 0.005)	0.948 (0.943 - 0.952)	0.000 (0.000 - 0.000)	0.996 (0.995 - 0.997)	0.000 (0.000 - 0.000)	-0.001 (-0.002 - -0.001)
	ECGCLIP-R18	0.003 (0.001 - 0.009)	0.969 (0.966 - 0.973)	0.000 (0.000 - 0.000)	0.997 (0.996 - 0.998)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - 0.000)
	ECGCLIP-R34	0.003 (0.001 - 0.010)	0.973 (0.967 - 0.978)	0.000 (0.000 - 0.000)	0.996 (0.995 - 0.997)	0.000 (0.000 - 0.000)	-0.001 (-0.002 - -0.001)
VPE	Merl-R18	0.003 (0.001 - 0.008)	0.967 (0.963 - 0.971)	0.000 (0.000 - 0.000)	0.995 (0.993 - 0.996)	0.000 (0.000 - 0.000)	-0.001 (-0.002 - -0.001)
	Random Init-R18	0.125 (0.000 - 1.000)	0.906 (0.807 - 1.000)	0.500 (0.000 - 1.000)	0.999 (0.999 - 1.000)	0.167 (0.000 - 0.500)	0.223 (-0.001 - 0.548)
	ECGCLIP-R18	0.125 (0.000 - 1.000)	0.873 (0.741 - 1.000)	0.500 (0.000 - 1.000)	0.999 (0.999 - 1.000)	0.167 (0.000 - 0.500)	0.223 (-0.001 - 0.548)
	ECGCLIP-R34	0.043 (0.001 - 0.250)	0.975 (0.948 - 1.000)	0.500 (0.000 - 1.000)	0.999 (0.999 - 1.000)	0.182 (0.000 - 0.500)	0.235 (0.000 - 0.577)
3° AVB	Merl-R18	0.050 (0.000 - 0.286)	0.907 (0.810 - 1.000)	0.500 (0.000 - 1.000)	0.999 (0.998 - 1.000)	0.154 (0.000 - 0.444)	0.213 (-0.001 - 0.522)
	Random Init-R18	0.036 (0.021 - 0.062)	0.727 (0.663 - 0.786)	0.024 (0.000 - 0.063)	0.997 (0.996 - 0.998)	0.035 (0.000 - 0.087)	0.034 (-0.005 - 0.092)
	ECGCLIP-R18	0.042 (0.023 - 0.075)	0.716 (0.653 - 0.775)	0.098 (0.037 - 0.167)	0.996 (0.995 - 0.997)	0.123 (0.048 - 0.202)	0.122 (0.045 - 0.204)
	ECGCLIP-R34	0.071 (0.036 - 0.137)	0.774 (0.714 - 0.827)	0.061 (0.013 - 0.118)	0.999 (0.998 - 0.999)	0.100 (0.022 - 0.184)	0.127 (0.029 - 0.230)
	Merl-R18	0.033 (0.020 - 0.055)	0.731 (0.670 - 0.786)	0.049 (0.011 - 0.100)	0.997 (0.995 - 0.998)	0.066 (0.016 - 0.131)	0.066 (0.011 - 0.136)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
2° 1 Type AVB	Random Init-R18	0.032 (0.020 - 0.054)	0.715 (0.662 - 0.763)	0.031 (0.000 - 0.069)	0.998 (0.997 - 0.998)	0.048 (0.000 - 0.105)	0.053 (-0.005 - 0.122)
	ECGCLIP-R18	0.037 (0.021 - 0.066)	0.644 (0.578 - 0.706)	0.083 (0.031 - 0.143)	0.995 (0.994 - 0.996)	0.105 (0.039 - 0.179)	0.103 (0.034 - 0.179)
	ECGCLIP-R34	0.046 (0.024 - 0.092)	0.708 (0.662 - 0.755)	0.073 (0.028 - 0.130)	0.996 (0.995 - 0.997)	0.098 (0.039 - 0.171)	0.098 (0.036 - 0.174)
2° 2 Type AVB	Merl-R18	0.030 (0.017 - 0.055)	0.660 (0.603 - 0.715)	0.031 (0.000 - 0.069)	0.996 (0.995 - 0.998)	0.044 (0.000 - 0.098)	0.043 (-0.006 - 0.103)
	Random Init-R18	0.017 (0.011 - 0.030)	0.590 (0.532 - 0.651)	0.010 (0.000 - 0.035)	0.999 (0.998 - 1.000)	0.019 (0.000 - 0.061)	0.029 (-0.004 - 0.099)
	ECGCLIP-R18	0.025 (0.014 - 0.057)	0.594 (0.534 - 0.664)	0.021 (0.000 - 0.054)	0.999 (0.998 - 0.999)	0.036 (0.000 - 0.090)	0.049 (-0.004 - 0.124)
LVH	ECGCLIP-R34	0.027 (0.011 - 0.075)	0.427 (0.355 - 0.502)	0.031 (0.000 - 0.074)	0.999 (0.999 - 1.000)	0.056 (0.000 - 0.126)	0.086 (-0.003 - 0.184)
	Merl-R18	0.015 (0.008 - 0.046)	0.461 (0.396 - 0.532)	0.010 (0.000 - 0.035)	0.999 (0.999 - 1.000)	0.019 (0.000 - 0.061)	0.029 (-0.004 - 0.100)
	Random Init-R18	0.471 (0.443 - 0.502)	0.784 (0.770 - 0.800)	0.315 (0.290 - 0.343)	0.979 (0.977 - 0.982)	0.430 (0.403 - 0.460)	0.417 (0.389 - 0.446)
LAH	ECGCLIP-R18	0.493 (0.464 - 0.523)	0.799 (0.785 - 0.814)	0.310 (0.284 - 0.336)	0.981 (0.979 - 0.984)	0.430 (0.400 - 0.458)	0.422 (0.393 - 0.450)
	ECGCLIP-R34	0.476 (0.447 - 0.507)	0.796 (0.783 - 0.811)	0.332 (0.308 - 0.359)	0.977 (0.974 - 0.980)	0.443 (0.416 - 0.471)	0.423 (0.394 - 0.452)
	Merl-R18	0.465 (0.437 - 0.497)	0.778 (0.764 - 0.794)	0.329 (0.305 - 0.356)	0.977 (0.974 - 0.980)	0.440 (0.414 - 0.469)	0.420 (0.392 - 0.452)
RVH	Random Init-R18	0.602 (0.567 - 0.638)	0.937 (0.930 - 0.945)	0.672 (0.641 - 0.700)	0.956 (0.952 - 0.960)	0.635 (0.610 - 0.658)	0.598 (0.571 - 0.623)
	ECGCLIP-R18	0.630 (0.596 - 0.667)	0.942 (0.934 - 0.949)	0.676 (0.647 - 0.707)	0.960 (0.956 - 0.964)	0.650 (0.626 - 0.674)	0.615 (0.589 - 0.641)
	ECGCLIP-R34	0.644 (0.609 - 0.679)	0.944 (0.937 - 0.951)	0.778 (0.751 - 0.805)	0.941 (0.936 - 0.946)	0.655 (0.632 - 0.678)	0.624 (0.600 - 0.649)
	Merl-R18	0.617 (0.581 - 0.652)	0.938 (0.929 - 0.945)	0.737 (0.710 - 0.762)	0.943 (0.939 - 0.948)	0.638 (0.615 - 0.661)	0.604 (0.578 - 0.628)
	Random Init-R18	0.043 (0.022 - 0.075)	0.723 (0.655 - 0.794)	0.211 (0.118 - 0.313)	0.986 (0.984 - 0.989)	0.133 (0.073 - 0.199)	0.135 (0.071 - 0.202)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
RAH	ECGCLIP-R18	0.044 (0.022 - 0.084)	0.724 (0.660 - 0.787)	0.169 (0.085 - 0.262)	0.990 (0.988 - 0.992)	0.130 (0.069 - 0.203)	0.126 (0.063 - 0.200)
	ECGCLIP-R34	0.047 (0.024 - 0.085)	0.737 (0.674 - 0.803)	0.197 (0.108 - 0.292)	0.988 (0.986 - 0.990)	0.136 (0.074 - 0.203)	0.135 (0.070 - 0.204)
	Merl-R18	0.047 (0.024 - 0.092)	0.725 (0.663 - 0.791)	0.169 (0.085 - 0.266)	0.992 (0.990 - 0.993)	0.143 (0.069 - 0.218)	0.138 (0.065 - 0.214)
	Random Init-R18	0.085 (0.055 - 0.127)	0.804 (0.764 - 0.840)	0.071 (0.032 - 0.116)	0.996 (0.995 - 0.997)	0.106 (0.049 - 0.167)	0.115 (0.049 - 0.181)
	ECGCLIP-R18	0.097 (0.063 - 0.147)	0.817 (0.781 - 0.850)	0.085 (0.041 - 0.131)	0.997 (0.996 - 0.998)	0.132 (0.066 - 0.197)	0.152 (0.076 - 0.227)
	ECGCLIP-R34	0.136 (0.091 - 0.202)	0.823 (0.785 - 0.859)	0.078 (0.037 - 0.128)	0.999 (0.998 - 1.000)	0.135 (0.067 - 0.211)	0.194 (0.107 - 0.280)
	Merl-R18	0.084 (0.051 - 0.129)	0.784 (0.743 - 0.823)	0.035 (0.007 - 0.072)	0.999 (0.998 - 1.000)	0.064 (0.014 - 0.126)	0.105 (0.024 - 0.191)
	Random Init-R18	0.136 (0.064 - 0.256)	0.854 (0.790 - 0.908)	0.327 (0.204 - 0.453)	0.991 (0.989 - 0.992)	0.206 (0.119 - 0.286)	0.216 (0.127 - 0.301)
	ECGCLIP-R18	0.251 (0.134 - 0.373)	0.880 (0.812 - 0.933)	0.385 (0.246 - 0.515)	0.994 (0.993 - 0.996)	0.305 (0.195 - 0.410)	0.308 (0.196 - 0.412)
	ECGCLIP-R34	0.431 (0.270 - 0.583)	0.898 (0.828 - 0.949)	0.615 (0.466 - 0.742)	0.995 (0.993 - 0.996)	0.467 (0.357 - 0.563)	0.478 (0.368 - 0.571)
AAI	Merl-R18	0.091 (0.046 - 0.177)	0.870 (0.810 - 0.919)	0.212 (0.105 - 0.328)	0.993 (0.992 - 0.995)	0.168 (0.085 - 0.258)	0.166 (0.081 - 0.257)
	Random Init-R18	0.270 (0.142 - 0.418)	0.924 (0.864 - 0.969)	0.356 (0.213 - 0.500)	0.997 (0.996 - 0.998)	0.368 (0.232 - 0.484)	0.365 (0.230 - 0.486)
	ECGCLIP-R18	0.300 (0.160 - 0.456)	0.917 (0.852 - 0.968)	0.289 (0.156 - 0.422)	0.998 (0.997 - 0.999)	0.347 (0.203 - 0.483)	0.352 (0.205 - 0.493)
	ECGCLIP-R34	0.219 (0.117 - 0.394)	0.908 (0.844 - 0.960)	0.133 (0.041 - 0.238)	0.999 (0.999 - 1.000)	0.211 (0.070 - 0.349)	0.257 (0.093 - 0.402)
VP	Merl-R18	0.302 (0.164 - 0.457)	0.923 (0.859 - 0.972)	0.356 (0.220 - 0.500)	0.997 (0.996 - 0.998)	0.352 (0.220 - 0.474)	0.349 (0.217 - 0.472)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S25: Detailed task-specific performance comparison on the external CPSC2018 cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
STEMI	Random Init-R18	0.111 (0.075 – 0.152)	0.647 (0.611 – 0.681)	0.055 (0.027 – 0.088)	0.999 (0.999 – 1.000)	0.102 (0.051 – 0.160)	0.197 (0.118 – 0.265)
	ECGCLIP-R18	0.123 (0.085 – 0.165)	0.681 (0.646 – 0.715)	0.055 (0.027 – 0.088)	0.999 (0.998 – 1.000)	0.100 (0.050 – 0.159)	0.179 (0.103 – 0.251)
	ECGCLIP-R34	0.153 (0.111 – 0.202)	0.739 (0.707 – 0.770)	0.068 (0.037 – 0.104)	0.998 (0.997 – 0.999)	0.120 (0.066 – 0.180)	0.179 (0.106 – 0.255)
	Merl-R18	0.101 (0.067 – 0.142)	0.598 (0.558 – 0.636)	0.050 (0.022 – 0.081)	0.998 (0.997 – 0.999)	0.091 (0.042 – 0.143)	0.147 (0.072 – 0.221)
NSR	Random Init-R18	0.358 (0.329 – 0.392)	0.778 (0.764 – 0.794)	0.948 (0.934 – 0.962)	0.371 (0.359 – 0.385)	0.314 (0.298 – 0.331)	0.231 (0.218 – 0.245)
	ECGCLIP-R18	0.319 (0.294 – 0.350)	0.774 (0.760 – 0.788)	0.942 (0.928 – 0.957)	0.382 (0.369 – 0.394)	0.316 (0.301 – 0.333)	0.233 (0.219 – 0.248)
	ECGCLIP-R34	0.301 (0.276 – 0.328)	0.769 (0.756 – 0.784)	0.955 (0.943 – 0.969)	0.371 (0.358 – 0.384)	0.316 (0.301 – 0.333)	0.236 (0.223 – 0.250)
	Merl-R18	0.372 (0.343 – 0.408)	0.783 (0.769 – 0.797)	0.939 (0.925 – 0.953)	0.381 (0.369 – 0.393)	0.315 (0.300 – 0.332)	0.230 (0.216 – 0.245)
VPB	Random Init-R18	0.640 (0.607 – 0.676)	0.846 (0.828 – 0.864)	0.593 (0.557 – 0.629)	0.968 (0.964 – 0.973)	0.633 (0.603 – 0.665)	0.596 (0.563 – 0.630)
	ECGCLIP-R18	0.677 (0.644 – 0.709)	0.857 (0.838 – 0.876)	0.641 (0.608 – 0.678)	0.970 (0.965 – 0.974)	0.672 (0.644 – 0.702)	0.638 (0.609 – 0.670)
	ECGCLIP-R34	0.731 (0.701 – 0.764)	0.883 (0.866 – 0.900)	0.743 (0.710 – 0.774)	0.967 (0.962 – 0.971)	0.729 (0.704 – 0.754)	0.698 (0.671 – 0.725)
	Merl-R18	0.647 (0.612 – 0.682)	0.854 (0.837 – 0.872)	0.594 (0.558 – 0.630)	0.970 (0.966 – 0.974)	0.640 (0.609 – 0.669)	0.604 (0.572 – 0.636)
APB	Random Init-R18	0.262 (0.232 – 0.297)	0.757 (0.739 – 0.775)	0.416 (0.377 – 0.454)	0.876 (0.867 – 0.883)	0.310 (0.282 – 0.340)	0.233 (0.200 – 0.265)
	ECGCLIP-R18	0.311 (0.277 – 0.352)	0.786 (0.768 – 0.803)	0.411 (0.369 – 0.449)	0.910 (0.903 – 0.917)	0.353 (0.323 – 0.383)	0.283 (0.250 – 0.315)
	ECGCLIP-R34	0.418 (0.377 – 0.465)	0.829 (0.812 – 0.846)	0.481 (0.440 – 0.519)	0.945 (0.940 – 0.951)	0.471 (0.438 – 0.504)	0.418 (0.383 – 0.453)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AF	Merl-R18	0.273 (0.242 – 0.309)	0.767 (0.750 – 0.785)	0.409 (0.368 – 0.446)	0.895 (0.887 – 0.901)	0.330 (0.301 – 0.358)	0.256 (0.223 – 0.287)
	Random Init-R18	0.875 (0.856 – 0.893)	0.949 (0.941 – 0.958)	0.536 (0.509 – 0.564)	0.995 (0.993 – 0.997)	0.687 (0.664 – 0.711)	0.678 (0.657 – 0.700)
	ECGCLIP-R18	0.867 (0.848 – 0.886)	0.941 (0.933 – 0.950)	0.515 (0.490 – 0.545)	0.997 (0.995 – 0.998)	0.673 (0.651 – 0.698)	0.669 (0.649 – 0.691)
	ECGCLIP-R34	0.851 (0.832 – 0.870)	0.937 (0.929 – 0.945)	0.518 (0.491 – 0.546)	0.997 (0.995 – 0.998)	0.676 (0.653 – 0.699)	0.672 (0.652 – 0.692)
RBBB	Merl-R18	0.875 (0.855 – 0.894)	0.946 (0.938 – 0.955)	0.559 (0.532 – 0.586)	0.996 (0.994 – 0.997)	0.707 (0.686 – 0.730)	0.697 (0.677 – 0.718)
	Random Init-R18	0.955 (0.946 – 0.964)	0.984 (0.981 – 0.987)	0.793 (0.774 – 0.811)	0.985 (0.981 – 0.988)	0.865 (0.853 – 0.877)	0.827 (0.812 – 0.841)
	ECGCLIP-R18	0.953 (0.941 – 0.963)	0.984 (0.981 – 0.987)	0.794 (0.774 – 0.812)	0.987 (0.983 – 0.990)	0.868 (0.855 – 0.879)	0.831 (0.816 – 0.845)
	ECGCLIP-R34	0.955 (0.945 – 0.964)	0.985 (0.982 – 0.987)	0.794 (0.775 – 0.811)	0.986 (0.982 – 0.989)	0.866 (0.854 – 0.877)	0.829 (0.814 – 0.842)
1° AVB	Merl-R18	0.954 (0.944 – 0.962)	0.983 (0.979 – 0.986)	0.786 (0.767 – 0.805)	0.985 (0.982 – 0.988)	0.861 (0.849 – 0.873)	0.822 (0.807 – 0.837)
	Random Init-R18	0.887 (0.864 – 0.908)	0.977 (0.970 – 0.982)	0.713 (0.679 – 0.747)	0.990 (0.987 – 0.992)	0.793 (0.768 – 0.816)	0.778 (0.751 – 0.802)
	ECGCLIP-R18	0.889 (0.865 – 0.909)	0.978 (0.972 – 0.984)	0.740 (0.708 – 0.771)	0.989 (0.986 – 0.992)	0.807 (0.783 – 0.831)	0.791 (0.765 – 0.816)
	ECGCLIP-R34	0.892 (0.868 – 0.914)	0.980 (0.975 – 0.985)	0.699 (0.664 – 0.732)	0.991 (0.989 – 0.993)	0.787 (0.762 – 0.811)	0.773 (0.747 – 0.797)
LBBB	Merl-R18	0.886 (0.862 – 0.908)	0.977 (0.970 – 0.983)	0.751 (0.719 – 0.781)	0.988 (0.985 – 0.991)	0.810 (0.787 – 0.832)	0.793 (0.768 – 0.816)
	Random Init-R18	0.903 (0.868 – 0.931)	0.979 (0.965 – 0.990)	0.627 (0.561 – 0.682)	0.999 (0.999 – 1.000)	0.763 (0.708 – 0.803)	0.776 (0.729 – 0.811)
	ECGCLIP-R18	0.908 (0.875 – 0.936)	0.983 (0.971 – 0.992)	0.606 (0.537 – 0.662)	1.000 (0.999 – 1.000)	0.749 (0.693 – 0.790)	0.765 (0.717 – 0.801)
	ECGCLIP-R34	0.899 (0.864 – 0.930)	0.980 (0.968 – 0.989)	0.627 (0.563 – 0.681)	0.999 (0.999 – 1.000)	0.763 (0.711 – 0.803)	0.776 (0.731 – 0.812)
	Merl-R18	0.909 (0.875 – 0.937)	0.983 (0.970 – 0.992)	0.610 (0.542 – 0.667)	1.000 (0.999 – 1.000)	0.752 (0.697 – 0.795)	0.767 (0.720 – 0.805)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S26: Detailed task-specific performance comparison on the external Shiyuan cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
NER	Random Init-R18	0.003 (0.003 - 0.003)	0.654 (0.635 - 0.675)	0.820 (0.783 - 0.853)	0.443 (0.441 - 0.445)	0.006 (0.005 - 0.007)	0.024 (0.020 - 0.027)
	ECGCLIP-R18	0.003 (0.003 - 0.003)	0.632 (0.612 - 0.654)	0.757 (0.720 - 0.796)	0.455 (0.453 - 0.457)	0.006 (0.005 - 0.006)	0.019 (0.015 - 0.023)
	ECGCLIP-R34	0.003 (0.002 - 0.003)	0.624 (0.603 - 0.646)	0.745 (0.706 - 0.785)	0.443 (0.441 - 0.445)	0.005 (0.005 - 0.006)	0.017 (0.013 - 0.021)
	Merl-R18	0.003 (0.002 - 0.003)	0.614 (0.593 - 0.635)	0.730 (0.688 - 0.771)	0.453 (0.451 - 0.455)	0.005 (0.005 - 0.006)	0.016 (0.013 - 0.020)
LAE	Random Init-R18	0.332 (0.326 - 0.337)	0.806 (0.803 - 0.808)	0.713 (0.708 - 0.719)	0.735 (0.733 - 0.737)	0.386 (0.382 - 0.389)	0.310 (0.306 - 0.314)
	ECGCLIP-R18	0.348 (0.342 - 0.354)	0.811 (0.809 - 0.813)	0.670 (0.665 - 0.676)	0.772 (0.770 - 0.774)	0.397 (0.393 - 0.401)	0.318 (0.313 - 0.322)
	ECGCLIP-R34	0.350 (0.344 - 0.356)	0.807 (0.805 - 0.809)	0.665 (0.660 - 0.671)	0.769 (0.767 - 0.771)	0.392 (0.388 - 0.396)	0.311 (0.307 - 0.316)
	Merl-R18	0.336 (0.331 - 0.342)	0.805 (0.803 - 0.807)	0.663 (0.657 - 0.669)	0.767 (0.765 - 0.769)	0.389 (0.386 - 0.393)	0.308 (0.304 - 0.312)
HCM	Random Init-R18	0.381 (0.373 - 0.389)	0.842 (0.839 - 0.845)	0.545 (0.538 - 0.552)	0.893 (0.892 - 0.895)	0.396 (0.391 - 0.401)	0.342 (0.337 - 0.348)
	ECGCLIP-R18	0.393 (0.386 - 0.401)	0.849 (0.846 - 0.852)	0.538 (0.530 - 0.545)	0.902 (0.901 - 0.903)	0.405 (0.400 - 0.411)	0.352 (0.347 - 0.358)
	ECGCLIP-R34	0.396 (0.388 - 0.404)	0.850 (0.848 - 0.853)	0.550 (0.543 - 0.557)	0.898 (0.896 - 0.899)	0.405 (0.400 - 0.411)	0.353 (0.347 - 0.359)
	Merl-R18	0.383 (0.376 - 0.391)	0.845 (0.843 - 0.848)	0.490 (0.483 - 0.497)	0.917 (0.916 - 0.918)	0.403 (0.398 - 0.409)	0.347 (0.341 - 0.354)
nHCM	Random Init-R18	0.348 (0.341 - 0.356)	0.840 (0.837 - 0.843)	0.560 (0.553 - 0.568)	0.883 (0.881 - 0.884)	0.382 (0.377 - 0.388)	0.332 (0.327 - 0.338)
	ECGCLIP-R18	0.357 (0.350 - 0.364)	0.846 (0.843 - 0.848)	0.578 (0.571 - 0.585)	0.881 (0.879 - 0.882)	0.388 (0.383 - 0.394)	0.341 (0.335 - 0.346)
	ECGCLIP-R34	0.356 (0.348 - 0.363)	0.847 (0.844 - 0.850)	0.608 (0.600 - 0.615)	0.866 (0.864 - 0.867)	0.382 (0.377 - 0.388)	0.338 (0.333 - 0.344)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
BAE	Merl-R18	0.348 (0.341 - 0.355)	0.843 (0.840 - 0.846)	0.530 (0.523 - 0.537)	0.899 (0.897 - 0.900)	0.390 (0.385 - 0.396)	0.338 (0.332 - 0.343)
	Random Init-R18	0.394 (0.378 - 0.411)	0.975 (0.973 - 0.977)	0.725 (0.710 - 0.739)	0.977 (0.976 - 0.978)	0.470 (0.459 - 0.482)	0.491 (0.480 - 0.502)
	ECGCLIP-R18	0.418 (0.401 - 0.435)	0.977 (0.976 - 0.979)	0.740 (0.724 - 0.754)	0.977 (0.977 - 0.978)	0.479 (0.467 - 0.491)	0.501 (0.489 - 0.512)
	ECGCLIP-R34	0.428 (0.412 - 0.445)	0.979 (0.977 - 0.980)	0.725 (0.710 - 0.739)	0.980 (0.979 - 0.980)	0.495 (0.483 - 0.507)	0.511 (0.500 - 0.523)
HFrfEF	Merl-R18	0.404 (0.387 - 0.421)	0.976 (0.974 - 0.977)	0.739 (0.725 - 0.755)	0.976 (0.975 - 0.976)	0.463 (0.452 - 0.474)	0.488 (0.477 - 0.499)
	Random Init-R18	0.489 (0.471 - 0.507)	0.949 (0.945 - 0.953)	0.625 (0.611 - 0.640)	0.985 (0.984 - 0.985)	0.510 (0.498 - 0.522)	0.508 (0.496 - 0.521)
	ECGCLIP-R18	0.570 (0.554 - 0.587)	0.957 (0.953 - 0.961)	0.552 (0.537 - 0.568)	0.992 (0.992 - 0.992)	0.557 (0.544 - 0.571)	0.549 (0.536 - 0.563)
	ECGCLIP-R34	0.581 (0.563 - 0.596)	0.962 (0.959 - 0.966)	0.612 (0.597 - 0.627)	0.990 (0.990 - 0.990)	0.569 (0.557 - 0.581)	0.562 (0.550 - 0.574)
RMWSSFLV	Merl-R18	0.549 (0.531 - 0.566)	0.954 (0.950 - 0.958)	0.596 (0.581 - 0.611)	0.990 (0.989 - 0.990)	0.554 (0.542 - 0.567)	0.547 (0.534 - 0.560)
	Random Init-R18	0.001 (0.000 - 0.006)	0.944 (0.847 - 0.998)	0.600 (0.000 - 1.000)	0.984 (0.984 - 0.985)	0.002 (0.000 - 0.004)	0.022 (0.000 - 0.040)
	ECGCLIP-R18	0.001 (0.000 - 0.004)	0.971 (0.940 - 0.997)	0.400 (0.000 - 1.000)	0.982 (0.981 - 0.982)	0.001 (0.000 - 0.002)	0.014 (-0.001 - 0.031)
	ECGCLIP-R34	0.001 (0.000 - 0.007)	0.961 (0.909 - 0.995)	0.400 (0.000 - 1.000)	0.981 (0.981 - 0.982)	0.001 (0.000 - 0.002)	0.013 (-0.001 - 0.030)
DCM	Merl-R18	0.001 (0.000 - 0.004)	0.969 (0.932 - 0.998)	0.400 (0.000 - 1.000)	0.984 (0.983 - 0.984)	0.001 (0.000 - 0.003)	0.014 (-0.001 - 0.033)
	Random Init-R18	0.002 (0.000 - 0.010)	0.992 (0.980 - 0.999)	0.667 (0.000 - 1.000)	0.994 (0.994 - 0.995)	0.003 (0.000 - 0.008)	0.032 (0.000 - 0.062)
	ECGCLIP-R18	0.005 (0.000 - 0.023)	0.997 (0.993 - 1.000)	1.000 (1.000 - 1.000)	0.992 (0.991 - 0.992)	0.003 (0.001 - 0.007)	0.040 (0.023 - 0.061)
	ECGCLIP-R34	0.011 (0.000 - 0.080)	0.996 (0.990 - 1.000)	0.667 (0.000 - 1.000)	0.993 (0.992 - 0.993)	0.003 (0.000 - 0.006)	0.029 (0.000 - 0.056)
	Merl-R18	0.006 (0.000 - 0.028)	0.995 (0.987 - 1.000)	0.667 (0.000 - 1.000)	0.995 (0.995 - 0.995)	0.004 (0.000 - 0.009)	0.034 (0.000 - 0.066)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
oHCM	Random Init-R18	0.033 (0.027 - 0.042)	0.837 (0.819 - 0.853)	0.165 (0.135 - 0.198)	0.995 (0.995 - 0.995)	0.105 (0.087 - 0.127)	0.109 (0.090 - 0.133)
	ECGCLIP-R18	0.035 (0.029 - 0.045)	0.833 (0.816 - 0.850)	0.150 (0.124 - 0.181)	0.996 (0.996 - 0.996)	0.109 (0.090 - 0.130)	0.110 (0.091 - 0.133)
	ECGCLIP-R34	0.037 (0.030 - 0.047)	0.846 (0.829 - 0.862)	0.207 (0.174 - 0.245)	0.994 (0.994 - 0.994)	0.118 (0.099 - 0.140)	0.127 (0.106 - 0.151)
	Merl-R18	0.034 (0.028 - 0.043)	0.837 (0.821 - 0.853)	0.181 (0.150 - 0.215)	0.995 (0.994 - 0.995)	0.111 (0.092 - 0.133)	0.117 (0.096 - 0.141)
LVDD	Random Init-R18	0.810 (0.808 - 0.813)	0.715 (0.713 - 0.717)	0.031 (0.030 - 0.032)	0.989 (0.988 - 0.989)	0.060 (0.059 - 0.062)	0.060 (0.057 - 0.063)
	ECGCLIP-R18	0.820 (0.818 - 0.823)	0.728 (0.726 - 0.731)	0.013 (0.012 - 0.013)	0.996 (0.996 - 0.997)	0.025 (0.024 - 0.026)	0.043 (0.040 - 0.046)
	ECGCLIP-R34	0.828 (0.826 - 0.831)	0.742 (0.740 - 0.744)	0.042 (0.041 - 0.043)	0.988 (0.987 - 0.988)	0.080 (0.078 - 0.082)	0.078 (0.075 - 0.081)
	Merl-R18	0.813 (0.811 - 0.815)	0.717 (0.715 - 0.719)	0.025 (0.025 - 0.026)	0.992 (0.991 - 0.993)	0.049 (0.048 - 0.051)	0.058 (0.055 - 0.061)
VA	Random Init-R18	0.343 (0.312 - 0.376)	0.980 (0.973 - 0.986)	0.444 (0.409 - 0.481)	0.998 (0.998 - 0.998)	0.444 (0.414 - 0.475)	0.442 (0.412 - 0.473)
	ECGCLIP-R18	0.351 (0.320 - 0.385)	0.983 (0.977 - 0.988)	0.429 (0.392 - 0.465)	0.998 (0.998 - 0.998)	0.424 (0.393 - 0.456)	0.422 (0.391 - 0.454)
	ECGCLIP-R34	0.380 (0.347 - 0.416)	0.985 (0.979 - 0.990)	0.388 (0.351 - 0.422)	0.999 (0.998 - 0.999)	0.426 (0.391 - 0.458)	0.427 (0.392 - 0.458)
	Merl-R18	0.324 (0.294 - 0.357)	0.981 (0.975 - 0.987)	0.492 (0.458 - 0.527)	0.997 (0.997 - 0.998)	0.438 (0.409 - 0.467)	0.438 (0.409 - 0.467)
RAE	Random Init-R18	0.018 (0.014 - 0.026)	0.949 (0.937 - 0.960)	0.069 (0.037 - 0.104)	0.997 (0.997 - 0.997)	0.035 (0.019 - 0.053)	0.038 (0.020 - 0.058)
	ECGCLIP-R18	0.022 (0.017 - 0.032)	0.956 (0.944 - 0.965)	0.468 (0.401 - 0.537)	0.980 (0.979 - 0.980)	0.043 (0.035 - 0.051)	0.099 (0.082 - 0.115)
	ECGCLIP-R34	0.023 (0.018 - 0.032)	0.953 (0.940 - 0.962)	0.440 (0.374 - 0.505)	0.981 (0.980 - 0.981)	0.042 (0.034 - 0.051)	0.095 (0.079 - 0.111)
	Merl-R18	0.021 (0.016 - 0.031)	0.954 (0.943 - 0.964)	0.106 (0.067 - 0.148)	0.996 (0.996 - 0.997)	0.044 (0.027 - 0.061)	0.052 (0.033 - 0.074)
AM	Random Init-R18	0.000 (0.000 - 0.000)	0.689 (0.583 - 0.781)	0.000 (0.000 - 0.000)	0.997 (0.997 - 0.997)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAM	ECGCLIP-R18	0.000 (0.000 - 0.001)	0.733 (0.639 - 0.816)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.001 (0.000 - 0.002)	0.712 (0.615 - 0.808)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Merl-R18	0.000 (0.000 - 0.000)	0.703 (0.605 - 0.793)	0.000 (0.000 - 0.000)	0.999 (0.998 - 0.999)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
	Random Init-R18	0.002 (0.001 - 0.008)	0.748 (0.680 - 0.813)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.009 (0.001 - 0.049)	0.731 (0.668 - 0.792)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	ECGCLIP-R34	0.004 (0.001 - 0.023)	0.753 (0.689 - 0.814)	0.055 (0.000 - 0.122)	0.999 (0.999 - 0.999)	0.025 (0.000 - 0.055)	0.029 (0.000 - 0.065)
	Merl-R18	0.004 (0.001 - 0.018)	0.772 (0.699 - 0.836)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.001 (0.000 - 0.001)	0.000 (0.000 - 0.000)
	Random Init-R18	0.248 (0.244 - 0.253)	0.756 (0.753 - 0.759)	0.679 (0.673 - 0.685)	0.693 (0.691 - 0.695)	0.322 (0.318 - 0.326)	0.242 (0.238 - 0.246)
AVC	ECGCLIP-R18	0.257 (0.252 - 0.262)	0.765 (0.762 - 0.767)	0.684 (0.678 - 0.690)	0.699 (0.697 - 0.701)	0.328 (0.324 - 0.331)	0.250 (0.246 - 0.255)
	ECGCLIP-R34	0.261 (0.256 - 0.266)	0.767 (0.764 - 0.770)	0.696 (0.690 - 0.701)	0.694 (0.692 - 0.696)	0.329 (0.325 - 0.333)	0.254 (0.249 - 0.258)
	Merl-R18	0.258 (0.252 - 0.262)	0.765 (0.762 - 0.768)	0.690 (0.685 - 0.696)	0.698 (0.696 - 0.700)	0.330 (0.326 - 0.333)	0.254 (0.249 - 0.258)
	Random Init-R18	0.403 (0.390 - 0.415)	0.904 (0.900 - 0.907)	0.585 (0.574 - 0.595)	0.961 (0.960 - 0.962)	0.433 (0.424 - 0.442)	0.424 (0.415 - 0.432)
TR	ECGCLIP-R18	0.442 (0.430 - 0.454)	0.911 (0.907 - 0.914)	0.577 (0.567 - 0.588)	0.967 (0.966 - 0.968)	0.460 (0.451 - 0.469)	0.447 (0.438 - 0.456)
	ECGCLIP-R34	0.445 (0.433 - 0.457)	0.913 (0.910 - 0.917)	0.587 (0.576 - 0.598)	0.966 (0.965 - 0.966)	0.458 (0.450 - 0.467)	0.447 (0.438 - 0.456)
	Merl-R18	0.427 (0.415 - 0.439)	0.905 (0.901 - 0.909)	0.541 (0.530 - 0.552)	0.971 (0.970 - 0.972)	0.457 (0.449 - 0.466)	0.441 (0.432 - 0.450)
MR	Random Init-R18	0.251 (0.239 - 0.263)	0.883 (0.878 - 0.887)	0.525 (0.511 - 0.538)	0.956 (0.955 - 0.957)	0.318 (0.309 - 0.327)	0.322 (0.313 - 0.331)
	ECGCLIP-R18	0.272 (0.260 - 0.285)	0.891 (0.886 - 0.896)	0.547 (0.535 - 0.560)	0.958 (0.957 - 0.958)	0.336 (0.327 - 0.345)	0.341 (0.332 - 0.350)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AR	ECGCLIP-R34	0.276 (0.264 - 0.288)	0.894 (0.889 - 0.898)	0.507 (0.494 - 0.521)	0.965 (0.965 - 0.966)	0.350 (0.340 - 0.359)	0.347 (0.337 - 0.357)
	Merl-R18	0.266 (0.254 - 0.278)	0.889 (0.884 - 0.893)	0.514 (0.501 - 0.527)	0.963 (0.962 - 0.964)	0.341 (0.332 - 0.350)	0.340 (0.331 - 0.350)
	Random Init-R18	0.118 (0.108 - 0.128)	0.827 (0.820 - 0.834)	0.163 (0.152 - 0.176)	0.989 (0.989 - 0.990)	0.180 (0.168 - 0.193)	0.169 (0.156 - 0.182)
	ECGCLIP-R18	0.146 (0.135 - 0.157)	0.841 (0.835 - 0.848)	0.182 (0.170 - 0.194)	0.990 (0.990 - 0.991)	0.206 (0.193 - 0.219)	0.196 (0.183 - 0.210)
	ECGCLIP-R34	0.148 (0.137 - 0.160)	0.842 (0.836 - 0.849)	0.159 (0.147 - 0.171)	0.992 (0.992 - 0.993)	0.196 (0.184 - 0.210)	0.192 (0.179 - 0.205)
	Merl-R18	0.139 (0.128 - 0.150)	0.838 (0.832 - 0.845)	0.171 (0.159 - 0.184)	0.990 (0.990 - 0.991)	0.193 (0.180 - 0.206)	0.183 (0.170 - 0.197)
MVPLAC	Random Init-R18	0.070 (0.065 - 0.075)	0.826 (0.821 - 0.832)	0.282 (0.268 - 0.299)	0.953 (0.952 - 0.954)	0.130 (0.123 - 0.138)	0.131 (0.123 - 0.140)
	ECGCLIP-R18	0.074 (0.069 - 0.079)	0.832 (0.826 - 0.837)	0.360 (0.344 - 0.377)	0.935 (0.934 - 0.936)	0.129 (0.123 - 0.136)	0.142 (0.134 - 0.150)
	ECGCLIP-R34	0.075 (0.070 - 0.080)	0.833 (0.828 - 0.839)	0.335 (0.320 - 0.352)	0.945 (0.945 - 0.946)	0.137 (0.130 - 0.144)	0.146 (0.137 - 0.155)
	Merl-R18	0.070 (0.066 - 0.076)	0.830 (0.825 - 0.836)	0.323 (0.309 - 0.339)	0.944 (0.943 - 0.945)	0.131 (0.124 - 0.137)	0.138 (0.130 - 0.146)
MS	Random Init-R18	0.133 (0.099 - 0.179)	0.962 (0.953 - 0.970)	0.366 (0.312 - 0.423)	0.997 (0.996 - 0.997)	0.185 (0.153 - 0.215)	0.211 (0.177 - 0.244)
	ECGCLIP-R18	0.182 (0.138 - 0.236)	0.966 (0.958 - 0.974)	0.376 (0.323 - 0.436)	0.997 (0.997 - 0.997)	0.210 (0.176 - 0.244)	0.233 (0.197 - 0.270)
	ECGCLIP-R34	0.198 (0.154 - 0.255)	0.967 (0.958 - 0.975)	0.352 (0.296 - 0.407)	0.998 (0.998 - 0.998)	0.239 (0.200 - 0.278)	0.251 (0.211 - 0.291)
AS	Merl-R18	0.155 (0.115 - 0.206)	0.964 (0.955 - 0.972)	0.369 (0.315 - 0.426)	0.997 (0.997 - 0.997)	0.199 (0.166 - 0.231)	0.222 (0.188 - 0.259)
	Random Init-R18	0.050 (0.038 - 0.064)	0.882 (0.868 - 0.896)	0.168 (0.138 - 0.201)	0.995 (0.995 - 0.996)	0.112 (0.091 - 0.134)	0.116 (0.095 - 0.138)
AS	ECGCLIP-R18	0.074 (0.055 - 0.095)	0.880 (0.866 - 0.896)	0.212 (0.180 - 0.247)	0.995 (0.995 - 0.996)	0.140 (0.118 - 0.162)	0.145 (0.124 - 0.168)
	ECGCLIP-R34	0.076 (0.058 - 0.099)	0.883 (0.868 - 0.898)	0.148 (0.119 - 0.180)	0.998 (0.998 - 0.998)	0.145 (0.119 - 0.175)	0.143 (0.116 - 0.173)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
BAV	Merl-R18	0.061 (0.045 - 0.079)	0.884 (0.870 - 0.898)	0.181 (0.151 - 0.212)	0.996 (0.995 - 0.996)	0.124 (0.104 - 0.146)	0.127 (0.106 - 0.151)
	Random Init-R18	0.005 (0.003 - 0.010)	0.651 (0.614 - 0.689)	0.056 (0.031 - 0.088)	0.997 (0.997 - 0.997)	0.031 (0.017 - 0.047)	0.033 (0.017 - 0.051)
	ECGCLIP-R18	0.006 (0.004 - 0.010)	0.670 (0.633 - 0.707)	0.072 (0.041 - 0.107)	0.997 (0.996 - 0.997)	0.036 (0.020 - 0.053)	0.040 (0.021 - 0.060)
	ECGCLIP-R34	0.006 (0.004 - 0.009)	0.684 (0.651 - 0.719)	0.068 (0.040 - 0.099)	0.997 (0.997 - 0.997)	0.036 (0.021 - 0.053)	0.039 (0.022 - 0.058)
PS	Merl-R18	0.006 (0.003 - 0.011)	0.648 (0.612 - 0.686)	0.068 (0.040 - 0.102)	0.996 (0.996 - 0.996)	0.029 (0.017 - 0.043)	0.034 (0.018 - 0.051)
	Random Init-R18	0.002 (0.001 - 0.007)	0.695 (0.599 - 0.796)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.004 (0.001 - 0.017)	0.802 (0.691 - 0.899)	0.065 (0.000 - 0.174)	0.999 (0.999 - 0.999)	0.015 (0.000 - 0.039)	0.023 (0.000 - 0.059)
	ECGCLIP-R34	0.007 (0.001 - 0.029)	0.743 (0.629 - 0.850)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
EA	Merl-R18	0.033 (0.000 - 0.083)	0.653 (0.544 - 0.765)	0.290 (0.140 - 0.464)	0.978 (0.978 - 0.979)	0.004 (0.002 - 0.006)	0.022 (0.009 - 0.035)
	Random Init-R18	0.000 (0.000 - 0.001)	0.847 (0.692 - 0.964)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.009 (0.001 - 0.046)	0.974 (0.943 - 0.996)	0.100 (0.000 - 0.333)	1.000 (1.000 - 1.000)	0.050 (0.000 - 0.158)	0.058 (0.000 - 0.188)
	ECGCLIP-R34	0.005 (0.001 - 0.020)	0.976 (0.939 - 0.995)	0.100 (0.000 - 0.333)	1.000 (1.000 - 1.000)	0.022 (0.000 - 0.072)	0.035 (0.000 - 0.115)
PR	Merl-R18	0.000 (0.000 - 0.001)	0.906 (0.824 - 0.966)	0.000 (0.000 - 0.000)	0.996 (0.996 - 0.997)	0.000 (0.000 - 0.000)	0.000 (-0.001 - 0.000)
	Random Init-R18	0.007 (0.004 - 0.011)	0.787 (0.746 - 0.827)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.005 (0.003 - 0.009)	0.794 (0.757 - 0.830)	0.025 (0.006 - 0.051)	1.000 (0.999 - 1.000)	0.030 (0.007 - 0.062)	0.030 (0.007 - 0.063)
	ECGCLIP-R34	0.008 (0.004 - 0.015)	0.809 (0.774 - 0.841)	0.018 (0.000 - 0.042)	1.000 (1.000 - 1.000)	0.029 (0.000 - 0.065)	0.034 (0.000 - 0.076)
	Merl-R18	0.016 (0.008 - 0.027)	0.871 (0.843 - 0.900)	0.294 (0.228 - 0.364)	0.986 (0.986 - 0.987)	0.030 (0.021 - 0.039)	0.065 (0.048 - 0.083)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
CHD	Random Init-R18	0.194 (0.164 - 0.227)	0.890 (0.876 - 0.904)	0.282 (0.251 - 0.312)	0.998 (0.998 - 0.998)	0.291 (0.263 - 0.321)	0.289 (0.261 - 0.319)
	ECGCLIP-R18	0.212 (0.182 - 0.245)	0.907 (0.896 - 0.919)	0.311 (0.280 - 0.346)	0.998 (0.997 - 0.998)	0.312 (0.283 - 0.342)	0.310 (0.281 - 0.340)
	ECGCLIP-R34	0.233 (0.200 - 0.268)	0.906 (0.893 - 0.918)	0.308 (0.276 - 0.341)	0.998 (0.998 - 0.998)	0.332 (0.302 - 0.365)	0.331 (0.301 - 0.364)
ASD	Merl-R18	0.196 (0.166 - 0.228)	0.898 (0.885 - 0.911)	0.292 (0.261 - 0.324)	0.997 (0.997 - 0.998)	0.286 (0.257 - 0.316)	0.284 (0.255 - 0.314)
	Random Init-R18	0.213 (0.179 - 0.253)	0.874 (0.857 - 0.890)	0.314 (0.278 - 0.351)	0.997 (0.997 - 0.998)	0.292 (0.262 - 0.325)	0.290 (0.260 - 0.323)
	ECGCLIP-R18	0.253 (0.216 - 0.292)	0.890 (0.875 - 0.905)	0.346 (0.310 - 0.382)	0.998 (0.998 - 0.998)	0.332 (0.300 - 0.363)	0.330 (0.298 - 0.362)
VSD	ECGCLIP-R34	0.258 (0.223 - 0.298)	0.891 (0.876 - 0.906)	0.390 (0.354 - 0.430)	0.997 (0.996 - 0.997)	0.305 (0.279 - 0.335)	0.310 (0.283 - 0.340)
	Merl-R18	0.222 (0.188 - 0.259)	0.885 (0.869 - 0.900)	0.394 (0.359 - 0.432)	0.996 (0.996 - 0.996)	0.291 (0.266 - 0.319)	0.299 (0.274 - 0.327)
	Random Init-R18	0.019 (0.010 - 0.042)	0.841 (0.810 - 0.869)	0.115 (0.067 - 0.170)	0.998 (0.998 - 0.999)	0.065 (0.036 - 0.094)	0.071 (0.040 - 0.104)
PFO	ECGCLIP-R18	0.028 (0.015 - 0.049)	0.875 (0.849 - 0.899)	0.216 (0.149 - 0.287)	0.997 (0.997 - 0.998)	0.079 (0.054 - 0.106)	0.101 (0.069 - 0.133)
	ECGCLIP-R34	0.042 (0.021 - 0.080)	0.890 (0.866 - 0.912)	0.180 (0.118 - 0.246)	0.998 (0.998 - 0.999)	0.096 (0.061 - 0.129)	0.107 (0.069 - 0.145)
	Merl-R18	0.024 (0.012 - 0.046)	0.852 (0.818 - 0.883)	0.151 (0.094 - 0.212)	0.997 (0.997 - 0.998)	0.058 (0.035 - 0.081)	0.072 (0.044 - 0.101)
PDA	Random Init-R18	0.006 (0.005 - 0.007)	0.604 (0.583 - 0.626)	0.042 (0.028 - 0.057)	0.992 (0.992 - 0.993)	0.024 (0.016 - 0.033)	0.022 (0.013 - 0.032)
	ECGCLIP-R18	0.006 (0.005 - 0.008)	0.627 (0.605 - 0.648)	0.050 (0.035 - 0.068)	0.990 (0.990 - 0.991)	0.024 (0.016 - 0.033)	0.023 (0.014 - 0.033)
	ECGCLIP-R34	0.007 (0.006 - 0.008)	0.638 (0.617 - 0.658)	0.067 (0.050 - 0.087)	0.987 (0.987 - 0.988)	0.027 (0.020 - 0.035)	0.028 (0.019 - 0.037)
	Merl-R18	0.006 (0.005 - 0.008)	0.617 (0.595 - 0.639)	0.057 (0.041 - 0.076)	0.991 (0.991 - 0.991)	0.030 (0.022 - 0.039)	0.029 (0.019 - 0.040)
	Random Init-R18	0.011 (0.003 - 0.037)	0.794 (0.723 - 0.854)	0.062 (0.013 - 0.127)	0.999 (0.999 - 1.000)	0.044 (0.010 - 0.090)	0.046 (0.011 - 0.093)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
TOF	ECGCLIP-R18	0.010 (0.003 - 0.027)	0.795 (0.726 - 0.858)	0.046 (0.000 - 0.108)	1.000 (1.000 - 1.000)	0.047 (0.000 - 0.105)	0.047 (0.000 - 0.105)
	ECGCLIP-R34	0.020 (0.006 - 0.061)	0.798 (0.732 - 0.859)	0.046 (0.000 - 0.108)	1.000 (1.000 - 1.000)	0.063 (0.000 - 0.137)	0.068 (0.000 - 0.149)
	Merl-R18	0.011 (0.003 - 0.032)	0.797 (0.733 - 0.857)	0.092 (0.032 - 0.175)	0.999 (0.999 - 0.999)	0.057 (0.019 - 0.105)	0.061 (0.020 - 0.112)
	Random Init-R18	0.037 (0.002 - 0.163)	0.868 (0.786 - 0.942)	0.370 (0.194 - 0.542)	0.992 (0.992 - 0.992)	0.011 (0.005 - 0.018)	0.045 (0.022 - 0.069)
	ECGCLIP-R18	0.010 (0.002 - 0.043)	0.856 (0.786 - 0.921)	0.556 (0.360 - 0.750)	0.925 (0.924 - 0.926)	0.002 (0.001 - 0.003)	0.020 (0.011 - 0.029)
	ECGCLIP-R34	0.013 (0.001 - 0.102)	0.872 (0.776 - 0.947)	0.704 (0.500 - 0.868)	0.918 (0.917 - 0.919)	0.002 (0.001 - 0.003)	0.025 (0.016 - 0.033)
	Merl-R18	0.002 (0.000 - 0.013)	0.766 (0.659 - 0.863)	0.111 (0.000 - 0.241)	0.998 (0.997 - 0.998)	0.011 (0.000 - 0.025)	0.025 (-0.001 - 0.054)
	Random Init-R18	0.392 (0.385 - 0.398)	0.784 (0.781 - 0.787)	0.505 (0.498 - 0.510)	0.876 (0.874 - 0.877)	0.420 (0.415 - 0.424)	0.330 (0.325 - 0.336)
PH	ECGCLIP-R18	0.399 (0.393 - 0.405)	0.790 (0.787 - 0.793)	0.462 (0.456 - 0.467)	0.902 (0.900 - 0.903)	0.425 (0.419 - 0.430)	0.339 (0.333 - 0.345)
	ECGCLIP-R34	0.404 (0.397 - 0.410)	0.791 (0.788 - 0.794)	0.451 (0.445 - 0.457)	0.908 (0.907 - 0.910)	0.427 (0.422 - 0.432)	0.344 (0.338 - 0.349)
	Merl-R18	0.394 (0.388 - 0.401)	0.786 (0.783 - 0.789)	0.494 (0.488 - 0.500)	0.883 (0.882 - 0.884)	0.422 (0.418 - 0.427)	0.334 (0.328 - 0.339)
	Random Init-R18	0.343 (0.337 - 0.348)	0.777 (0.774 - 0.779)	0.346 (0.341 - 0.351)	0.917 (0.916 - 0.919)	0.362 (0.357 - 0.367)	0.274 (0.268 - 0.280)
AD	ECGCLIP-R18	0.373 (0.368 - 0.379)	0.798 (0.795 - 0.801)	0.390 (0.385 - 0.396)	0.911 (0.910 - 0.913)	0.391 (0.386 - 0.397)	0.303 (0.297 - 0.308)
	ECGCLIP-R34	0.381 (0.374 - 0.387)	0.801 (0.799 - 0.804)	0.457 (0.451 - 0.463)	0.888 (0.887 - 0.890)	0.412 (0.407 - 0.417)	0.318 (0.312 - 0.323)
	Merl-R18	0.364 (0.358 - 0.371)	0.791 (0.789 - 0.794)	0.402 (0.397 - 0.408)	0.903 (0.902 - 0.905)	0.390 (0.385 - 0.395)	0.298 (0.292 - 0.304)
Dex	Random Init-R18	0.001 (0.001 - 0.001)	0.740 (0.682 - 0.799)	0.000 (0.000 - 0.000)	0.998 (0.998 - 0.998)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)
	ECGCLIP-R18	0.085 (0.045 - 0.140)	0.958 (0.938 - 0.977)	0.538 (0.433 - 0.656)	0.994 (0.993 - 0.994)	0.055 (0.040 - 0.072)	0.123 (0.096 - 0.154)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
PE	ECGCLIP-R34	0.159 (0.082 – 0.257)	0.956 (0.937 – 0.973)	0.359 (0.253 – 0.468)	0.999 (0.999 – 0.999)	0.186 (0.125 – 0.248)	0.212 (0.145 – 0.276)
	Merl-R18	0.008 (0.004 – 0.016)	0.876 (0.830 – 0.917)	0.282 (0.186 – 0.385)	0.992 (0.991 – 0.992)	0.023 (0.015 – 0.034)	0.057 (0.037 – 0.079)
	Random Init-R18	0.209 (0.202 – 0.217)	0.755 (0.750 – 0.760)	0.241 (0.233 – 0.249)	0.967 (0.966 – 0.968)	0.265 (0.257 – 0.273)	0.229 (0.221 – 0.237)
	ECGCLIP-R18	0.228 (0.220 – 0.236)	0.766 (0.761 – 0.770)	0.279 (0.271 – 0.288)	0.963 (0.962 – 0.963)	0.288 (0.280 – 0.296)	0.249 (0.241 – 0.258)
	ECGCLIP-R34	0.230 (0.222 – 0.237)	0.767 (0.762 – 0.771)	0.265 (0.256 – 0.273)	0.966 (0.966 – 0.967)	0.286 (0.277 – 0.294)	0.249 (0.240 – 0.257)
	Merl-R18	0.216 (0.209 – 0.224)	0.759 (0.754 – 0.763)	0.267 (0.258 – 0.275)	0.962 (0.961 – 0.963)	0.275 (0.267 – 0.283)	0.235 (0.227 – 0.244)
CP	Random Init-R18	0.003 (0.000 – 0.023)	0.965 (0.886 – 1.000)	0.250 (0.000 – 1.000)	0.998 (0.998 – 0.998)	0.004 (0.000 – 0.013)	0.022 (0.000 – 0.066)
	ECGCLIP-R18	0.003 (0.000 – 0.012)	0.991 (0.982 – 0.999)	0.000 (0.000 – 0.000)	1.000 (1.000 – 1.000)	0.000 (0.000 – 0.000)	0.000 (0.000 – 0.000)
	ECGCLIP-R34	0.006 (0.000 – 0.031)	0.993 (0.979 – 1.000)	0.250 (0.000 – 0.800)	1.000 (0.999 – 1.000)	0.020 (0.000 – 0.066)	0.051 (0.000 – 0.146)
	Merl-R18	0.003 (0.000 – 0.012)	0.987 (0.967 – 0.999)	0.250 (0.000 – 1.000)	0.999 (0.999 – 0.999)	0.011 (0.000 – 0.036)	0.038 (0.000 – 0.111)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S27: Detailed task-specific performance comparison on the external BIDMC cohort.

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
LAE	Random Init-R18	0.683 (0.679 – 0.687)	0.697 (0.695 – 0.700)	0.341 (0.338 – 0.345)	0.848 (0.845 – 0.851)	0.461 (0.458 – 0.465)	0.218 (0.213 – 0.223)
	ECGCLIP-R18	0.687 (0.683 – 0.691)	0.700 (0.698 – 0.703)	0.264 (0.261 – 0.268)	0.891 (0.888 – 0.893)	0.388 (0.384 – 0.392)	0.197 (0.192 – 0.202)
	ECGCLIP-R34	0.695 (0.691 – 0.699)	0.709 (0.706 – 0.712)	0.222 (0.219 – 0.225)	0.914 (0.912 – 0.917)	0.342 (0.338 – 0.346)	0.187 (0.183 – 0.192)
	Merl-R18	0.689 (0.685 – 0.693)	0.702 (0.699 – 0.704)	0.260 (0.257 – 0.263)	0.892 (0.890 – 0.895)	0.383 (0.380 – 0.387)	0.195 (0.190 – 0.200)
BAE	Random Init-R18	0.599 (0.593 – 0.604)	0.711 (0.708 – 0.715)	0.215 (0.211 – 0.219)	0.962 (0.961 – 0.964)	0.335 (0.330 – 0.340)	0.282 (0.276 – 0.287)
	ECGCLIP-R18	0.615 (0.610 – 0.620)	0.726 (0.722 – 0.729)	0.226 (0.222 – 0.230)	0.963 (0.962 – 0.965)	0.350 (0.345 – 0.355)	0.297 (0.292 – 0.302)
	ECGCLIP-R34	0.622 (0.617 – 0.627)	0.734 (0.731 – 0.737)	0.184 (0.180 – 0.187)	0.974 (0.973 – 0.975)	0.298 (0.293 – 0.304)	0.275 (0.270 – 0.280)
	Merl-R18	0.610 (0.605 – 0.615)	0.718 (0.715 – 0.722)	0.225 (0.221 – 0.229)	0.963 (0.962 – 0.964)	0.348 (0.343 – 0.353)	0.295 (0.289 – 0.300)
HFrfEF	Random Init-R18	0.443 (0.422 – 0.465)	0.854 (0.845 – 0.862)	0.473 (0.453 – 0.493)	0.945 (0.942 – 0.948)	0.485 (0.469 – 0.503)	0.428 (0.410 – 0.448)
	ECGCLIP-R18	0.511 (0.488 – 0.530)	0.857 (0.848 – 0.865)	0.419 (0.397 – 0.439)	0.969 (0.967 – 0.972)	0.497 (0.476 – 0.516)	0.461 (0.439 – 0.480)
	ECGCLIP-R34	0.532 (0.509 – 0.553)	0.875 (0.866 – 0.883)	0.440 (0.419 – 0.458)	0.968 (0.966 – 0.970)	0.512 (0.493 – 0.530)	0.474 (0.454 – 0.492)
	Merl-R18	0.486 (0.464 – 0.506)	0.852 (0.843 – 0.861)	0.461 (0.440 – 0.480)	0.958 (0.955 – 0.960)	0.504 (0.485 – 0.521)	0.455 (0.435 – 0.474)
oHCM	Random Init-R18	0.073 (0.067 – 0.080)	0.664 (0.653 – 0.673)	0.049 (0.042 – 0.057)	0.993 (0.992 – 0.993)	0.076 (0.065 – 0.087)	0.076 (0.063 – 0.090)
	ECGCLIP-R18	0.064 (0.059 – 0.070)	0.653 (0.642 – 0.663)	0.030 (0.025 – 0.037)	0.996 (0.996 – 0.996)	0.052 (0.043 – 0.062)	0.063 (0.051 – 0.078)
	ECGCLIP-R34	0.072 (0.066 – 0.079)	0.654 (0.643 – 0.664)	0.057 (0.049 – 0.065)	0.993 (0.992 – 0.993)	0.087 (0.076 – 0.098)	0.088 (0.075 – 0.101)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
VA	Merl-R18	0.063 (0.058 - 0.068)	0.656 (0.646 - 0.666)	0.040 (0.033 - 0.047)	0.994 (0.993 - 0.994)	0.064 (0.053 - 0.075)	0.065 (0.052 - 0.078)
	Random Init-R18	0.948 (0.927 - 0.967)	0.714 (0.657 - 0.766)	0.047 (0.032 - 0.063)	0.988 (0.957 - 1.000)	0.090 (0.062 - 0.118)	0.052 (-0.001 - 0.082)
	ECGCLIP-R18	0.951 (0.929 - 0.968)	0.718 (0.660 - 0.772)	0.054 (0.038 - 0.072)	0.988 (0.959 - 1.000)	0.103 (0.073 - 0.134)	0.059 (0.013 - 0.088)
	ECGCLIP-R34	0.946 (0.922 - 0.965)	0.704 (0.644 - 0.761)	0.022 (0.011 - 0.034)	0.988 (0.958 - 1.000)	0.043 (0.023 - 0.065)	0.020 (-0.052 - 0.058)
RAE	Merl-R18	0.944 (0.920 - 0.962)	0.691 (0.628 - 0.748)	0.069 (0.050 - 0.090)	0.975 (0.937 - 1.000)	0.129 (0.096 - 0.165)	0.055 (0.003 - 0.095)
	Random Init-R18	0.581 (0.576 - 0.586)	0.622 (0.618 - 0.625)	0.021 (0.020 - 0.022)	0.994 (0.993 - 0.995)	0.041 (0.038 - 0.043)	0.066 (0.061 - 0.072)
	ECGCLIP-R18	0.603 (0.598 - 0.608)	0.641 (0.638 - 0.645)	0.129 (0.126 - 0.132)	0.966 (0.965 - 0.967)	0.221 (0.216 - 0.225)	0.178 (0.173 - 0.184)
	ECGCLIP-R34	0.609 (0.604 - 0.614)	0.643 (0.640 - 0.647)	0.137 (0.133 - 0.140)	0.967 (0.966 - 0.968)	0.232 (0.227 - 0.237)	0.191 (0.185 - 0.196)
TR	Merl-R18	0.606 (0.601 - 0.611)	0.638 (0.634 - 0.641)	0.016 (0.015 - 0.017)	0.997 (0.996 - 0.997)	0.032 (0.030 - 0.034)	0.069 (0.063 - 0.074)
	Random Init-R18	0.343 (0.336 - 0.351)	0.798 (0.794 - 0.801)	0.648 (0.641 - 0.656)	0.790 (0.788 - 0.793)	0.397 (0.392 - 0.403)	0.319 (0.313 - 0.325)
	ECGCLIP-R18	0.375 (0.367 - 0.383)	0.809 (0.806 - 0.813)	0.614 (0.606 - 0.621)	0.827 (0.825 - 0.829)	0.417 (0.411 - 0.422)	0.337 (0.331 - 0.344)
	ECGCLIP-R34	0.376 (0.368 - 0.384)	0.812 (0.809 - 0.816)	0.617 (0.609 - 0.624)	0.830 (0.828 - 0.832)	0.422 (0.416 - 0.427)	0.343 (0.337 - 0.350)
MR	Merl-R18	0.362 (0.354 - 0.370)	0.802 (0.799 - 0.806)	0.583 (0.575 - 0.590)	0.840 (0.838 - 0.842)	0.414 (0.408 - 0.420)	0.332 (0.325 - 0.338)
	Random Init-R18	0.272 (0.266 - 0.279)	0.762 (0.758 - 0.766)	0.483 (0.475 - 0.491)	0.842 (0.840 - 0.844)	0.352 (0.346 - 0.358)	0.259 (0.252 - 0.266)
	ECGCLIP-R18	0.288 (0.281 - 0.295)	0.772 (0.769 - 0.776)	0.469 (0.461 - 0.477)	0.859 (0.857 - 0.861)	0.362 (0.356 - 0.368)	0.270 (0.263 - 0.277)
	ECGCLIP-R34	0.296 (0.290 - 0.304)	0.777 (0.773 - 0.781)	0.449 (0.441 - 0.457)	0.872 (0.871 - 0.874)	0.365 (0.358 - 0.371)	0.274 (0.266 - 0.281)
	Merl-R18	0.284 (0.277 - 0.290)	0.768 (0.764 - 0.772)	0.467 (0.459 - 0.476)	0.856 (0.854 - 0.858)	0.358 (0.352 - 0.364)	0.265 (0.259 - 0.272)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
AR	Random Init-R18	0.107 (0.100 - 0.115)	0.645 (0.634 - 0.655)	0.111 (0.100 - 0.122)	0.969 (0.967 - 0.970)	0.135 (0.122 - 0.148)	0.098 (0.085 - 0.112)
	ECGCLIP-R18	0.111 (0.104 - 0.120)	0.651 (0.641 - 0.661)	0.124 (0.113 - 0.135)	0.966 (0.964 - 0.967)	0.145 (0.133 - 0.157)	0.105 (0.092 - 0.119)
	ECGCLIP-R34	0.115 (0.107 - 0.124)	0.658 (0.648 - 0.669)	0.107 (0.097 - 0.118)	0.974 (0.972 - 0.975)	0.138 (0.125 - 0.151)	0.107 (0.094 - 0.121)
	Merl-R18	0.110 (0.103 - 0.118)	0.646 (0.636 - 0.656)	0.135 (0.124 - 0.147)	0.958 (0.956 - 0.960)	0.146 (0.134 - 0.159)	0.101 (0.088 - 0.113)
MS	Random Init-R18	0.111 (0.088 - 0.138)	0.753 (0.729 - 0.776)	0.166 (0.135 - 0.200)	0.988 (0.986 - 0.989)	0.193 (0.158 - 0.228)	0.180 (0.144 - 0.217)
	ECGCLIP-R18	0.127 (0.100 - 0.158)	0.749 (0.726 - 0.772)	0.172 (0.141 - 0.206)	0.989 (0.988 - 0.991)	0.208 (0.173 - 0.244)	0.199 (0.162 - 0.236)
	ECGCLIP-R34	0.130 (0.105 - 0.160)	0.763 (0.741 - 0.785)	0.150 (0.120 - 0.183)	0.993 (0.992 - 0.994)	0.207 (0.167 - 0.245)	0.212 (0.172 - 0.251)
	Merl-R18	0.117 (0.093 - 0.149)	0.739 (0.715 - 0.762)	0.182 (0.149 - 0.215)	0.987 (0.986 - 0.989)	0.206 (0.170 - 0.241)	0.193 (0.156 - 0.229)
AS	Random Init-R18	0.527 (0.514 - 0.541)	0.606 (0.597 - 0.614)	0.051 (0.046 - 0.056)	0.977 (0.974 - 0.980)	0.094 (0.085 - 0.104)	0.077 (0.061 - 0.092)
	ECGCLIP-R18	0.528 (0.516 - 0.541)	0.608 (0.600 - 0.617)	0.070 (0.065 - 0.076)	0.967 (0.963 - 0.970)	0.126 (0.117 - 0.136)	0.085 (0.070 - 0.100)
	ECGCLIP-R34	0.536 (0.524 - 0.549)	0.614 (0.605 - 0.622)	0.047 (0.042 - 0.051)	0.982 (0.979 - 0.984)	0.087 (0.079 - 0.096)	0.082 (0.066 - 0.096)
	Merl-R18	0.524 (0.512 - 0.537)	0.606 (0.598 - 0.615)	0.075 (0.069 - 0.081)	0.966 (0.962 - 0.970)	0.133 (0.123 - 0.144)	0.091 (0.076 - 0.107)
PS	Random Init-R18	0.003 (0.002 - 0.005)	0.631 (0.567 - 0.687)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	ECGCLIP-R18	0.002 (0.001 - 0.002)	0.433 (0.371 - 0.497)	0.000 (0.000 - 0.000)	0.999 (0.999 - 1.000)	0.000 (0.000 - 0.000)	-0.001 (-0.001 - -0.001)
	ECGCLIP-R34	0.002 (0.001 - 0.003)	0.530 (0.468 - 0.590)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)
	Merl-R18	0.002 (0.002 - 0.003)	0.536 (0.464 - 0.602)	0.043 (0.000 - 0.092)	0.976 (0.974 - 0.977)	0.006 (0.000 - 0.014)	0.005 (-0.007 - 0.019)
PR	Random Init-R18	0.036 (0.033 - 0.040)	0.621 (0.609 - 0.634)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	0.000 (0.000 - 0.000)

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
VSD	ECGCLIP-R18	0.021 (0.020 - 0.022)	0.481 (0.467 - 0.494)	0.002 (0.001 - 0.005)	0.998 (0.998 - 0.998)	0.004 (0.001 - 0.009)	0.001 (-0.005 - 0.008)
	ECGCLIP-R34	0.015 (0.014 - 0.016)	0.396 (0.383 - 0.408)	0.000 (0.000 - 0.000)	1.000 (1.000 - 1.000)	0.000 (0.000 - 0.000)	-0.002 (-0.002 - -0.001)
	Merl-R18	0.041 (0.038 - 0.046)	0.675 (0.664 - 0.687)	0.087 (0.076 - 0.099)	0.968 (0.967 - 0.969)	0.065 (0.057 - 0.074)	0.043 (0.034 - 0.052)
	Random Init-R18	0.797 (0.735 - 0.857)	0.613 (0.550 - 0.687)	0.091 (0.055 - 0.126)	0.945 (0.891 - 0.989)	0.163 (0.104 - 0.221)	0.057 (-0.053 - 0.138)
	ECGCLIP-R18	0.774 (0.708 - 0.841)	0.581 (0.508 - 0.657)	0.142 (0.097 - 0.185)	0.890 (0.824 - 0.950)	0.240 (0.170 - 0.301)	0.041 (-0.068 - 0.142)
	ECGCLIP-R34	0.793 (0.728 - 0.850)	0.584 (0.513 - 0.660)	0.154 (0.108 - 0.199)	0.934 (0.880 - 0.978)	0.261 (0.191 - 0.326)	0.115 (0.019 - 0.192)
	Merl-R18	0.784 (0.720 - 0.848)	0.589 (0.519 - 0.663)	0.091 (0.057 - 0.128)	0.945 (0.892 - 0.988)	0.163 (0.106 - 0.221)	0.057 (-0.050 - 0.138)
	Random Init-R18	0.311 (0.305 - 0.318)	0.751 (0.747 - 0.755)	0.798 (0.793 - 0.804)	0.580 (0.577 - 0.584)	0.376 (0.371 - 0.380)	0.268 (0.263 - 0.273)
	ECGCLIP-R18	0.330 (0.323 - 0.337)	0.764 (0.761 - 0.768)	0.745 (0.738 - 0.752)	0.652 (0.649 - 0.656)	0.395 (0.390 - 0.400)	0.286 (0.281 - 0.292)
	ECGCLIP-R34	0.333 (0.326 - 0.340)	0.765 (0.762 - 0.769)	0.693 (0.686 - 0.700)	0.697 (0.694 - 0.700)	0.401 (0.395 - 0.405)	0.288 (0.282 - 0.293)
PH	Merl-R18	0.317 (0.310 - 0.324)	0.753 (0.749 - 0.757)	0.792 (0.786 - 0.799)	0.587 (0.584 - 0.590)	0.377 (0.372 - 0.381)	0.269 (0.263 - 0.274)
	Random Init-R18	0.286 (0.281 - 0.292)	0.678 (0.674 - 0.682)	0.155 (0.149 - 0.160)	0.945 (0.944 - 0.947)	0.217 (0.210 - 0.223)	0.145 (0.138 - 0.152)
	ECGCLIP-R18	0.302 (0.296 - 0.308)	0.692 (0.688 - 0.696)	0.192 (0.187 - 0.198)	0.931 (0.930 - 0.933)	0.250 (0.243 - 0.256)	0.161 (0.154 - 0.168)
	ECGCLIP-R34	0.309 (0.303 - 0.315)	0.696 (0.692 - 0.699)	0.249 (0.243 - 0.255)	0.910 (0.908 - 0.911)	0.292 (0.286 - 0.299)	0.184 (0.177 - 0.191)
AD	Merl-R18	0.296 (0.291 - 0.302)	0.693 (0.689 - 0.696)	0.228 (0.222 - 0.234)	0.915 (0.913 - 0.916)	0.275 (0.270 - 0.282)	0.170 (0.164 - 0.178)
	Random Init-R18	0.146 (0.132 - 0.160)	0.718 (0.705 - 0.729)	0.840 (0.823 - 0.858)	0.446 (0.440 - 0.452)	0.156 (0.149 - 0.163)	0.136 (0.126 - 0.145)
	ECGCLIP-R18	0.145 (0.132 - 0.159)	0.721 (0.709 - 0.732)	0.822 (0.802 - 0.840)	0.477 (0.471 - 0.484)	0.161 (0.153 - 0.168)	0.141 (0.131 - 0.151)
PE							

Disease	Model	PRAUC	ROAUC	Sensitivity	Specificity	F1 Score	MCC
	ECGCLIP-R34	0.140 (0.129 - 0.151)	0.723 (0.710 - 0.734)	0.819 (0.801 - 0.837)	0.487 (0.481 - 0.493)	0.162 (0.155 - 0.170)	0.144 (0.134 - 0.153)
	Merl-R18	0.142 (0.130 - 0.156)	0.714 (0.701 - 0.725)	0.815 (0.795 - 0.833)	0.468 (0.461 - 0.473)	0.157 (0.149 - 0.164)	0.133 (0.123 - 0.142)

Metrics include both threshold-independent measures (PRAUC, ROAUC) and threshold-dependent measures (Sensitivity, Specificity, F1 Score, MCC). Data are presented as point estimates followed by 95% CIs in parentheses.

Table S28: Data efficiency analysis comparing PRAUC performance across varying training data ratios.

Cohort	Training	Random Init-R18	ECGCLIP-R18	ECGCLIP-R34	Merl-R18
	Data Ratio				
Zhongshan	1%	0.186±0.255	<u>0.244±0.284</u>	0.289±0.303	0.212±0.268
	10%	0.369±0.321	<u>0.403±0.318</u>	0.436±0.320	0.381±0.321
	100%	0.392±0.325	<u>0.433±0.324</u>	0.482±0.326	0.404±0.326
Xiamen	1%	0.232±0.301	<u>0.275±0.316</u>	0.320±0.334	0.254±0.314
	10%	0.394±0.341	<u>0.432±0.336</u>	0.462±0.329	0.400±0.343
	100%	0.402±0.349	<u>0.446±0.340</u>	0.491±0.338	0.405±0.346
MIMIC-IV-ECG	1%	0.173±0.244	<u>0.189±0.253</u>	0.208±0.269	0.164±0.227
	10%	<u>0.251±0.302</u>	0.249±0.299	0.255±0.302	0.241±0.295
	100%	0.259±0.309	<u>0.264±0.310</u>	0.273±0.315	0.248±0.299
UKB	1%	0.213±0.316	<u>0.245±0.332</u>	0.266±0.334	0.225±0.322
	10%	0.318±0.350	<u>0.333±0.354</u>	0.344±0.360	0.322±0.357
	100%	0.332±0.354	<u>0.354±0.362</u>	0.367±0.366	0.332±0.356
Chapman	1%	0.169±0.250	<u>0.188±0.269</u>	0.212±0.290	0.173±0.253
	10%	0.255±0.293	<u>0.258±0.293</u>	0.260±0.293	0.254±0.292
	100%	0.267±0.302	<u>0.275±0.303</u>	0.290±0.304	0.266±0.301
PTBXL	1%	0.247±0.255	<u>0.277±0.281</u>	0.300±0.290	0.264±0.268
	10%	0.405±0.321	0.389±0.319	0.389±0.318	<u>0.395±0.332</u>
	100%	0.407±0.323	0.405±0.326	<u>0.406±0.324</u>	0.399±0.334
Georgia	1%	0.175±0.247	<u>0.198±0.267</u>	0.225±0.282	0.189±0.261
	10%	0.270±0.300	<u>0.282±0.300</u>	0.284±0.302	0.272±0.300
	100%	0.282±0.307	<u>0.296±0.310</u>	0.310±0.315	0.280±0.311
CPSC2018	1%	0.525±0.330	<u>0.547±0.335</u>	0.568±0.338	0.541±0.335

Cohort	Training				
	Data Ratio	Random Init-R18	ECGCLIP-R18	ECGCLIP-R34	Merl-R18
	10%	<u>0.617±0.336</u>	0.610±0.336	0.629±0.317	0.611±0.342
	100%	0.624±0.335	<u>0.631±0.330</u>	0.650±0.313	0.627±0.334
	1%	0.240	<u>0.270</u>	0.298	0.253
Average	10%	0.360	<u>0.370</u>	0.382	0.360
	100%	0.370	<u>0.388</u>	0.409	0.370

Reported values represent the mean \pm standard deviation, averaged across all tasks for each cohort. **Bold** denotes the highest PRAUC score, and underlines denotes the second-highest score within each data ratio setting.