

ECG-LLM: Foundation Model for ECG-Based Cardiac Reasoning

Alexander Selivanov^{*1}, Friederike Jungmann^{1,3}, Jan Kehrer², Karl-Ludwig Laugwitz²,
Eimo Martens², and Daniel Rueckert^{1,4,5}

¹Chair for AI in Healthcare and Medicine, Technical University of Munich (TUM) and
TUM University Hospital, Germany

²Department of Cardiology, TUM University Hospital, Germany

³Institute for Diagnostic and Interventional Radiology, TUM University Hospital,
Technical University of Munich (TUM), Munich, Germany

⁴Department of Computing, Imperial College London, UK

⁵Munich Center for Machine Learning (MCML), Germany

Abstract

Electrocardiography (ECG) is an inexpensive, standard-of-care test for cardiac symptoms, but front-line triage often lacks immediate access to definitive imaging such as echocardiography (ECHO) or cardiac magnetic resonance (CMR). Furthermore, most existing ECG-AI systems are limited to fixed diagnostic labels or automated reports, constraining their use for patient-specific clinical reasoning. To address this gap, we introduce ECG-LLM, an ECG-conditioned large language model trained across four cohorts comprising 679,112 ECG studies from 186,409 patients. Using a novel multimodal-to-language supervision strategy, ECG-LLM is trained on clinically structured question-answer pairs derived from ECG signals, clinical context, CMR, and ECHO. This unified approach enables the model to answer diverse cardiovascular questions from a 12-lead ECG alone, spanning both conventional interpretation and phenotypes not directly visible on standard ECGs. ECG-LLM successfully recovers conventional ECG measurements, such as heart rate, and strongly predicts complex CMR-derived phenotypes, including ventricular and atrial volumes and ventricular function. Crucially, it detects vital echocardiographic phenotypes, including increased LV wall thickness, aortic stenosis, and right-ventricular systolic dysfunction. On standard ECG understanding tasks, ECG-LLM matches or exceeds existing baselines for diagnostic report generation and the ECG-QA benchmark. By moving beyond fixed-label prediction, this multimodal framework provides clinically valuable, question-driven cardiovascular reasoning to support general practitioner and front-line triage decisions when specialist review is delayed.

1 Introduction

Cardiovascular diseases (CVDs) remain the leading cause of death worldwide [1, 2], yet early detection in front-line care and general practitioner (GP) practice is often constrained by limited access to advanced diagnostic tools such as imaging. In primary care settings, triage decisions frequently rely on clinical history and a 12-lead electrocardiogram (ECG) as the first objective cardiac test [3]. While imaging modalities such as echocardiography (ECHO) and cardiac magnetic resonance (CMR) provide the "gold standard" quantitative assessment of cardiac structure and function [4, 5, 6, 7, 8, 9], they require specialized equipment, trained operators, and expert interpretation. They are also often limited by high costs, long waiting lists, and, in the case of CMR, longer workflows and patient contraindications. This creates a practical clinical

^{*}alexander.selivanov@tum.de

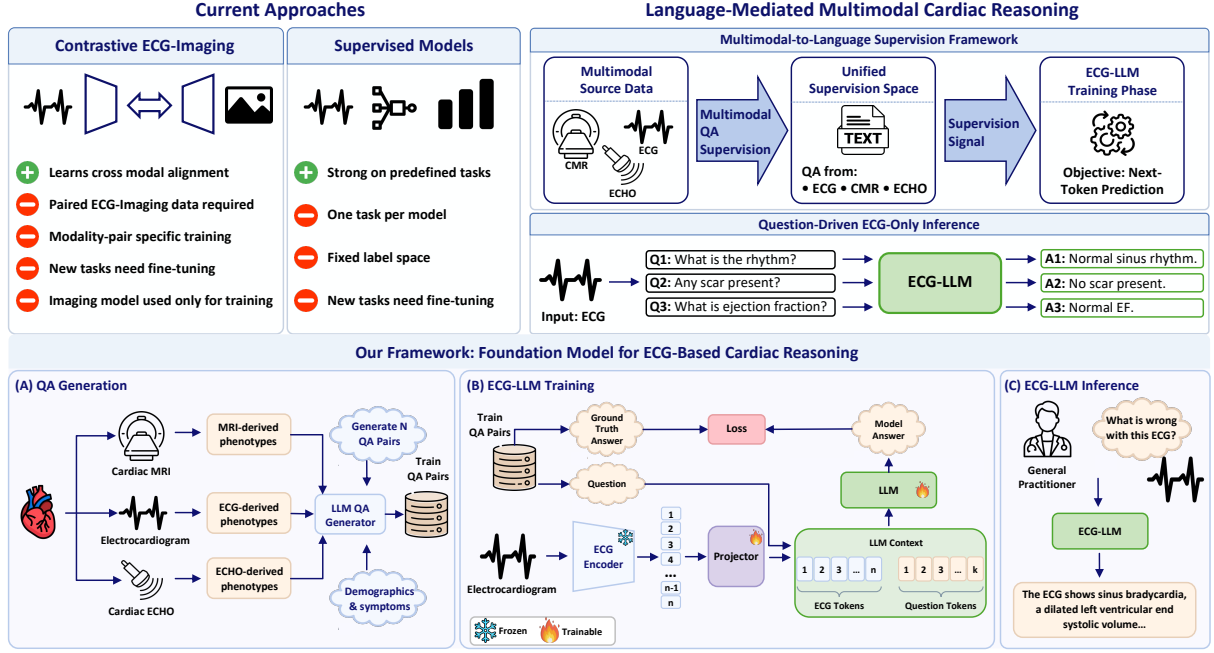


Figure 1: Overview of language-mediated multimodal cardiac reasoning. Upper panels: existing ECG AI approaches and our proposed multimodal QA supervision framework. Lower panels: (A) generation of QA pairs from ECG, CMR, ECHO, and clinical metadata; (B) ECG-LLM training with ECG tokens and question–answer supervision; (C) ECG-only inference for free-text clinical reasoning.

limitation: first-line triage often relies on the immediately available ECG, while the imaging-grade information needed for clinical decisions may be delayed until scheduled echocardiography or CMR is performed.

In this work, we address this gap with ECG-LLM, an ECG-conditioned large language model (LLM) that operates directly on the raw multichannel 12-lead ECG waveform, rather than on a rendered ECG image, and enables free-text cardiac question answering from the ECG alone. Our central hypothesis is that ECG waveforms, echocardiography, CMR, and clinical labels/reports are different observations of a shared underlying cardiovascular state. We therefore propose to use natural language as a unified supervision space for representing this state. Natural language provides a dense and flexible supervision format: it can describe clinically meaningful findings that rely on multiple tests, make heterogeneous information comparable through question–answer (QA) pairs, and match the natural way clinicians ask questions and communicate diagnoses. Instead of training separate models to predict a limited set of predefined categories, as in standard supervised approaches, we use language-based QA supervision to train a single ECG-conditioned LLM to answer clinical questions about both ECG findings and broader imaging-related cardiac phenotypes. Figure 1 summarizes this transition from conventional ECG AI toward language-mediated multimodal cardiac reasoning. The upper panels contrast existing AI approaches with our multimodal QA supervision framework, while the lower panels show how multimodal cardiac information is converted into QA supervision, used to train ECG-LLM, and applied for ECG-only clinical reasoning at inference time.

We build this framework across four large cohorts: UK Biobank [10], MIMIC-IV-ECG [11], PTB-XL [12], and EchoNext [13]. For each ECG study, we collect structured metadata describing all available ECG measurements and labels, imaging-derived phenotypes, patient demographics, symptoms, and clinical context when available. We then use this metadata to generate over 5 million clinically grounded QA pairs, formulated as natural GP-style queries about ECG findings and broader cardiac phenotypes. Using these QA pairs, we train ECG-

LLM, an autoregressive large language model conditioned on a 12-lead ECG raw signal and available patient context, to answer cardiac questions in free text. At inference time, ECG-LLM operates solely on the 12-lead ECG and patient context, without echocardiography, CMR, or task-specific classification heads, supporting ECG-based triage and second-opinion use in settings where imaging is delayed or unavailable. We describe the overall framework in more detail in Section 4.1, the construction of clinical QA supervision from cardiac records in Section 4.4, and we show representative examples of ECG-LLM answering diverse questions about ECG findings and CMR- and ECHO-derived cardiac phenotypes in Section 2.6.

The motivation for ECG-LLM is grounded in a broader shift in ECG AI: deep learning models are no longer limited to detecting visible waveform abnormalities, but can also infer clinically relevant cardiac phenotypes that are not directly measured on the ECG. Deep learning has substantially improved automated ECG interpretation, evolving from direct waveform analysis toward the inference of hidden phenotypes. Early models reached expert-level performance in arrhythmia detection and multi-label diagnosis [14, 15, 16]. Building on this, supervised approaches demonstrated that the 12-lead ECG can detect "indirect" structural diseases—such as aortic stenosis, reduced LVEF, and hypertrophic cardiomyopathy, as well as future arrhythmia risk [17, 18, 19, 20, 21, 22]. Further studies indicate that ECG representations can recover diastolic function, support unsupervised disease profiling, and act as digital biomarkers for hypertension-related risk [23, 24, 25].

Despite this progress, most supervised ECG models are designed to answer only the specific diagnostic questions defined before training. For example, a model trained to detect reduced left ventricle ejection fraction (LVEF) or aortic stenosis usually returns only that predefined output, rather than answering a clinician’s broader question in natural language. This makes such models difficult to use for open-ended triage, where the relevant question may vary across patients and clinical settings. Multimodal contrastive methods partially address this limitation by aligning ECG and imaging representations to transfer broader cardiovascular knowledge across modalities [26, 27]. However, these methods usually learn representations rather than directly answering clinical questions. In practice, they still require additional training of a prediction model to extract each specific phenotype of interest. This limits their ability to support direct, question-driven inference at the point of care. Their scalability is also constrained by reliance on imaging data, which remains costly to acquire, difficult to standardize, and often restricted from sharing due to privacy and infrastructure limitations.

In parallel, multimodal vision–language models (VLMs) have moved medical AI toward conversational reasoning and free-text generation [28, 29, 30, 31, 32]. Recent ECG-focused vision–language models shift toward language-based interpretation [33, 34, 35, 36], representing an important step toward clinically usable ECG-language systems. In practice, however, these models are primarily trained to reproduce ECG reports or generate textual interpretations of waveform-level findings, tasks that largely mirror existing automated reporting systems installed in ECG machines. As a result, their outputs remain focused on ECG-derived observations and do not directly target the broader cardiac phenotypes typically assessed using imaging, nor the targeted clinical questions encountered during triage. We build on this progress by extending ECG-language supervision beyond standard ECG machine-derived output: ECG-LLM uses multimodal cardiac information as dense QA supervision so that an ECG-conditioned LLM can answer flexible, question-driven queries about both ECG findings and imaging-related phenotypes from the raw ECG signal alone. The goal is to support triage and provide a second-opinion tool in settings where imaging is unavailable, delayed, or time- or access-constrained.

In summary, our contributions are:

- A multimodal-to-language supervision strategy that shifts ECG-language training beyond reproducing ECG-machine outputs toward question-driven reasoning over both ECG-derived and imaging-related cardiac phenotypes.

- A large-scale 12-lead ECG instruction-tuning resource and QA-generation pipeline, built from unified ECG measurements, diagnostic reports, ECHO and CMR phenotypes, clinical context, and 5M+ clinically realistic GP-style QA pairs. We plan to release the QA-generation pipeline and all shareable derived resources, subject to cohort-specific data-use agreements.
- ECG-LLM, an open-source ECG-conditioned autoregressive model that directly processes raw 12-lead ECG waveforms and answers free-text cardiac questions from ECG alone, without rendered ECG images, imaging at inference, or task-specific classification heads.
- A systematic evaluation across UK Biobank, MIMIC-IV-ECG, PTB-XL, and EchoNext, covering ECG interpretation, CMR- and ECHO-derived phenotype inference, the ECG-QA benchmark, and diagnostic report generation.

2 Results

2.1 Evaluation framework

ECG-LLM is an autoregressive large language model (LLM) conditioned on 12-lead ECGs and trained with multimodal question–answer supervision, as described in Methods (Section 4.1, Figure 1). At inference, the model receives a 5-s centre crop of the 12-lead ECG recording and a natural-language question, then generates a free-text answer (Figure 1C). We evaluate this ECG-to-text LLM without retrieval-augmented generation, auxiliary classification heads, or direct access to imaging measurements. Predictions are therefore conditioned only on the learned ECG representation and the question. We intentionally exclude retrieval-augmented generation because it can shift ECG interpretation from the ECG-conditioned language model to an external retrieval or prediction module. This setup tests how much clinically relevant information can be represented and queried through a generative ECG-conditioned language model. We evaluate five capabilities: (i) ECG-related question answering, (ii) CMR-related phenotype question answering from ECG only, (iii) echocardiography-related phenotype question answering from ECG only, (iv) ECG report generation, and (v) standard benchmark performance on ECG-QA [37].

2.2 ECG-LLM answers ECG and CMR phenotype questions from ECG alone

We evaluated whether ECG-LLM can answer clinical questions about both conventional ECG-derived phenotypes and CMR-derived imaging phenotypes using the ECG as the only input. This setting reflects a first-line care use case in which a clinician may ask targeted questions about rhythm, conduction, repolarisation, intervals, and imaging-derived structural or functional phenotypes, such as chamber size, ventricular function, myocardial wall thickness, strain, and possible remodelling, before deciding whether further imaging is needed. Existing ECG-language models such as ECG-Chat [33], PULSE-ECG [34], and ECG-Doctor [35] represent important steps toward ECG-language interpretation, but they primarily focus on diagnostic interpretation or report generation. Here, the model is instead queried in free text about both ECG measurements and imaging-derived cardiac phenotypes. We evaluated this task on 1,000 UK Biobank [10] test ECGs sampled from unique patients with non-trivial ECG device diagnoses, defined as any diagnosis other than the normal pair "Normal sinus rhythm" and "Normal ECG". This sampling avoids an evaluation dominated by completely normal ECGs, although class imbalance remains.

We organized ECG-derived and CMR-derived targets into clinically connected question groups, each with paraphrased natural-language questions and prespecified phenotypes. For

ECG-LLM Performance on UK Biobank ECG- and CMR-derived phenotypes

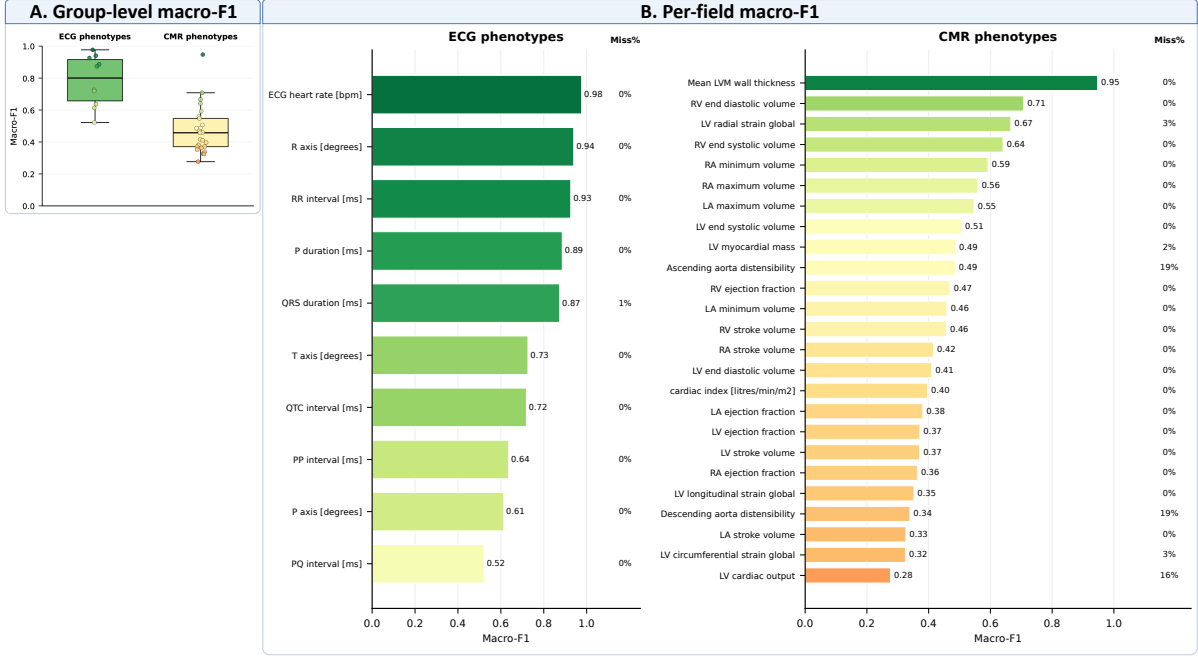


Figure 2: ECG-only question-answering performance for ECG-derived and CMR-derived phenotype categories in UK Biobank. The figure shows group-level and per-field class-macro F1 for ECG and CMR-derived phenotypes. Miss % indicates the fraction of cases in which no valid category was extracted from the generated answer for that phenotype.

example, a heart-rate group can target both heart rate and RR interval. For each ECG and question group, we sampled one question, generated an answer, and extracted only the phenotypes assigned to that group. This kept the evaluation aligned with the clinical intent of the question and avoided rewarding or penalizing unrelated statements. The ECG-derived phenotypes include heart rate, RR/PP/PQ intervals, P duration, QRS duration, QTc interval, and frontal P/R/T axes. The CMR-derived fields, derived from short-axis and long-axis CMR views, included ventricular and atrial volumes, ejection fractions, stroke volumes, LV myocardial mass, global wall thickness, strain, cardiac output, cardiac index, and aortic distensibility. The complete list of evaluated fields is provided in Supplementary Table 13, and the question groups, questions per group, and target phenotypes are listed in Supplementary Table 7. Because ECG-LLM generates free-text answers, we mapped each answer to the predefined category space using a separate schema-constrained extractor LLM. Details and the prompt used are provided in Supplementary Section 3. The extracted category for each phenotype was compared with the ground-truth category, treating each field as a single-label multi-class classification task. We report class-macro F1 as the primary metric because the category distributions are strongly imbalanced, with many fields dominated by normal values. Accuracy, balanced accuracy, weighted F1, coverage, missing rate, majority-class baseline performance, number of classes, and full class distributions are reported in Supplementary Table 8.

Figure 2 shows a higher class-macro F1 for ECG-derived fields than for CMR-derived fields. The strongest ECG-derived results are observed for directly ECG-visible properties, including heart rate, RR interval, R axis, P duration, and QRS duration. T axis and QTc interval show intermediate performance, while P axis and PP/PQ intervals are weaker, suggesting that atrial timing and atrial orientation are less consistently recovered than rate-, axis-, and ventricular conduction-related fields. Several CMR-derived phenotypes nevertheless show a clear signal above the majority baseline, especially LV mean myocardial wall thickness, RV end-diastolic volume, LV radial strain, RV end-systolic volume, RA minimum and maximum volume, LA

Confusion Matrices for UK Biobank ECG- and CMR-derived phenotypes

ECG-derived phenotypes									
ECG heart rate [bpm]					R axis [degrees]				
marked bradycardia	150	7	0	0	0	markedly short	2	0	0
bradycardia	0	443	17	0	0	short	0	2	0
normal	1	17	361	0	0	normal	0	1	369
mild tachycardia	0	0	0	2	0	prolonged	0	0	13
marked tachycardia	0	0	0	0	2	markedly prolonged	0	0	0
P duration [ms]					QRS duration [ms]				
normal	908	0			normal	896	12	1	
					mildly prolonged	8	22	0	
prolonged	25	46			prolonged	2	1	45	
T axis [degrees]									
leftward	14	16	12	2					
borderline	2	124	49	1					
normal	4	12	718	2					
rightward	2	4	6	32					
CMR-derived phenotypes									
Mean LVM wall thickness					RV end diastolic volume				
normal	762	33			small	1	0	7	
thinned	2	193			normal	0	959	0	
					dilated	0	0	32	
RV end systolic volume					RA minimum volume				
small	0	0	5		small	10	6	0	
normal	0	964	0		normal	7	874	0	
dilated	0	0	30		dilated	4	72	15	
					RA maximum volume				
					small	64	25	0	
					normal	38	830	0	
					dilated	0	30	1	

Figure 3: Confusion matrices for the top six ECG-derived and top six CMR-derived phenotype fields ranked by macro-F1 in the UK Biobank free-form phenotype extraction evaluation. Rows show ground-truth categories and columns show predicted categories.

maximum volume, and LV end-systolic volume. The strong result for LV wall thickness may partly reflect the training setup, because ECG-LLM was trained on QA pairs from ECG, CMR, and ECHO phenotypes, and wall-thickness-related remodelling appears in both the CMR and ECHO phenotype tasks. In contrast, weaker performance is observed for cardiac output, cardiac index, aortic distensibility, longitudinal and circumferential strain, and stroke-volume-related fields. These variables reflect hemodynamic function, vascular properties, or derived mechanical measurements rather than features directly measured on the ECG. For some phenotypes, performance remains close to the majority baseline, suggesting a weak ECG signal or that the model may partially collapse toward the most common category, leveraging its prior textual knowledge.

The confusion matrices in Figure 3A–B show errors for top-performing ECG-derived and CMR-derived fields. For ECG-derived variables, the strongest fields show clear diagonal structure, especially for heart rate, RR interval, R axis, P duration, and QRS duration. Errors are mainly concentrated in neighbouring or clinically related categories, although T-axis categories

remain less stable. For CMR-derived variables, the matrices are more heterogeneous. Mean LVM wall thickness shows strong separation between normal and thinned categories, while ventricular volume fields show useful separation of dominant normal or dilated categories. However, rare classes are less stable and are sometimes shifted toward more common categories. This partly reflects the difficulty of forcing continuous phenotypes, particularly values close to category boundaries, into strict discrete labels. In such borderline cases, the model may assign an incorrect adjacent category. However, because the answer is generated in free text, the model can still mention that the finding is borderline or that the prediction is uncertain. Thus, the CMR results support the presence of coarse structural signal in the ECG representation, while also showing that rare categories and precise severity boundaries remain challenging.

Taken together, this evaluation shows that ECG-LLM can answer free-text clinical questions about both conventional ECG measurements and CMR-derived imaging phenotypes from the ECG representation. The strongest performance is observed for ECG-proximal measurements, confirming that the model has learned clinically meaningful ECG concepts, while the above-baseline CMR results show that it also captures coarse structural signals that are not part of standard ECG-machine output. This distinguishes ECG-LLM from report-centred ECG-language systems: it is not limited to describing the ECG itself, but can describe imaging-related cardiac phenotypes that may help determine whether further cardiac imaging is warranted.

2.3 ECG-only question answering for echocardiographic phenotypes

To test whether ECG-LLM can reason about echocardiographic phenotypes from ECG alone, we evaluated it on the public EchoNext ECG–echocardiogram dataset [13]. The evaluation used 5,442 paired ECG–ECHO studies from the EchoNext test set. For each study, the model received only the 12-lead ECG and a phenotype-specific natural-language question. Because ECG-LLM was trained to produce free-text clinical answers rather than fixed binary predictions, we used a fixed two-turn protocol: first, the model answered the clinical question in free text. Second, within the same conversation, it was asked to provide its final answer as "yes" or "no". This preserved the model’s natural response format while producing a standardized binary label for evaluation. The final answer was parsed and compared with the echocardiographic ground truth. The exact questions used for each phenotype are provided in the Supplementary Table 9. We evaluated eleven ECHO-derived binary labels covering clinically thresholded or report-graded abnormalities (Supplementary Table 13): reduced left ventricular ejection fraction ($LVEF \leq 45\%$), increased left ventricular wall thickness ($LVWT \geq 1.3$ cm), moderate-or-greater right-ventricular systolic dysfunction, moderate-or-greater aortic stenosis, aortic regurgitation, mitral regurgitation, tricuspid regurgitation, pulmonary regurgitation, moderate-or-large pericardial effusion, elevated pulmonary artery systolic pressure ($PASP \geq 45$ mmHg), and elevated tricuspid regurgitation maximum velocity ($TR V_{max} \geq 3.2$ m/s).

Table 1: ECG-LLM performance on EchoNext binary echocardiographic phenotypes. All tasks were evaluated on 5,442 ECG–echocardiogram pairs using a fixed ECG-only yes/no question-answering protocol. Prev. denotes positive-class prevalence; Pos. rate denotes the fraction of ECG-LLM predictions parsed as positive.

Phenotype	Prev.	Pos. rate	Acc.	Prec.	Rec.	F1
LVWT ≥ 1.3 cm	19.5	18.5	0.911	0.787	0.748	0.767
AS $\geq$ moderate	5.3	8.2	0.956	0.552	0.857	0.671
RV systolic dysf. $\geq$ moderate	7.7	7.5	0.940	0.615	0.599	0.607
PASP ≥ 45 mmHg	12.8	9.2	0.908	0.697	0.498	0.581
LVEF $\leq 45\%$	17.7	32.6	0.774	0.424	0.783	0.550
TR Vmax ≥ 3.2 m/s	6.9	13.5	0.875	0.293	0.573	0.388
MR $\geq$ moderate	6.2	20.1	0.834	0.242	0.786	0.370
TR $\geq$ moderate	6.5	23.0	0.813	0.236	0.836	0.368
AR $\geq$ moderate	1.2	0.9	0.985	0.367	0.273	0.313
Pericardial eff. $\geq$ moderate/large	1.3	4.5	0.948	0.061	0.217	0.095
PR $\geq$ moderate	0.4	0.5	0.992	0.037	0.050	0.043

As shown in Table 1, ECG-LLM recovered a clear ECG-derived signal for several structural and functional echocardiographic phenotypes. The strongest result was observed for increased LV wall thickness, with an F1-score of 0.77, precision of 0.79, and recall of 0.75. The model also performed well for moderate-or-greater aortic stenosis, reaching an F1-score of 0.67 and a recall of 0.86 despite a positive prevalence of only 5.3%. Ventricular systolic dysfunction was also captured: moderate-or-greater right-ventricular systolic dysfunction reached an F1-score of 0.61 with balanced precision and recall, while reduced LVEF reached an F1-score of 0.55 with high recall of 0.78. Pulmonary-pressure-related phenotypes showed moderate performance. PASP ≥ 45 mmHg reached an F1-score of 0.58, with a precision of 0.70 and a recall of 0.50, whereas TR Vmax ≥ 3.2 m/s reached an F1-score of 0.39. Regurgitant valve phenotypes were more difficult. Moderate-or-greater mitral and tricuspid regurgitation had high recall, 0.79 and 0.84, respectively, but low precision around 0.24, resulting in F1-scores of 0.37. Aortic regurgitation reached an F1-score of 0.31. Thus, for several valve-related phenotypes, ECG-LLM often identified abnormal cases but also overcalled positives. The weakest results were observed for pulmonary regurgitation and pericardial effusion. These phenotypes were rare in the test set, with prevalences of 0.4% and 1.3%, and achieved F1-scores of 0.04 and 0.10, respectively. Accuracy was high for these tasks, but this was driven mainly by the dominant negative class. We also trained a supervised linear probe baseline using the same frozen ECG encoder as in ECG-LLM. On the EchoNext test set, this supervised ViT achieved a macro-AUROC of 0.79, macro-AUPRC of 0.24, and macro-F1 of 0.27, confirming that the encoder itself captures the echocardiographic signal. At the phenotype level, it performed best on reduced LVEF, increased LV wall thickness, and right-ventricular systolic dysfunction, consistent with the strongest ECG-LLM phenotypes.

Phenotype	Columbia Mini-Model [13]			Supervised ViT			ECG-LLM		
	AUROC	AUPRC	F1	AUROC	AUPRC	F1	P	R	F1
LVWT ≥ 1.3 cm	0.73	0.37	0.45	0.75	0.42	0.47	0.79	0.75	0.77
AS $\geq$ moderate	0.86	0.25	0.25	0.79	0.19	0.21	0.55	0.86	0.67
RV dysf $\geq$ moderate	0.87	0.43	0.37	0.86	0.43	0.38	0.62	0.60	0.61
PASP ≥ 45 mmHg	0.77	0.36	0.39	0.73	0.29	0.36	0.70	0.50	0.58
LVEF $\leq 45\%$	0.85	0.60	0.55	0.87	0.62	0.60	0.42	0.78	0.55
TR Vmax ≥ 3.2 m/s	0.75	0.21	0.24	0.71	0.15	0.23	0.29	0.57	0.39
MR $\geq$ moderate	0.81	0.23	0.26	0.80	0.24	0.27	0.24	0.79	0.37
TR $\geq$ moderate	0.83	0.30	0.29	0.80	0.23	0.28	0.24	0.84	0.37
AR $\geq$ moderate	0.74	0.04	0.05	0.74	0.03	0.05	0.37	0.27	0.31
Pericardial eff. $\geq$ mod	0.77	0.08	0.05	0.75	0.04	0.06	0.06	0.22	0.10
PR $\geq$ moderate	0.83	0.15	0.03	0.88	0.03	0.03	0.04	0.05	0.04

Table 2: Per-phenotype comparison on the EchoNext public test set. Columbia Mini-Model results are taken from the original EchoNext study [13]. Supervised ViT denotes a linear probe trained on top of the same frozen ECG encoder used in ECG-LLM. ECG-LLM predicts binary labels through ECG-only free-text question answering followed by a yes/no final-answer prompt and post hoc parsing. For Columbia Mini-Model and Supervised ViT, we report AUROC, AUPRC, and F1. For ECG-LLM, we report precision, recall, and F1 from parsed yes/no outputs. Bold indicates the highest F1-score for each phenotype.

The EchoNext benchmark was released together with the Columbia Mini-Model, whose performance on the same public test set was reported in the original study [13]. As summarized in Table 2, ECG-LLM achieved the highest F1-score for 10 of 11 phenotypes among the evaluated approaches, including LV wall thickness, aortic stenosis, right-ventricular systolic dysfunction, PASP elevation, TR Vmax elevation, mitral regurgitation, tricuspid regurgitation, aortic regurgitation, pericardial effusion, and pulmonary regurgitation. The only phenotype for which the supervised ViT achieved a higher F1-score was reduced LVEF, reaching 0.60 compared with 0.55 for ECG-LLM. Together, these results show that ECG-only echocardiographic question answering captures clinically meaningful signals for selected phenotypes, especially LV wall thickness, aortic stenosis, right-ventricular systolic dysfunction, PASP elevation, and reduced LVEF. Performance was weaker but still non-random for mitral and tricuspid regurgitation, and remained unreliable for rare or weakly ECG-expressed findings such as pulmonary regurgitation and pericardial effusion. These findings strengthen ECG-LLM as an ECG-based screening and reasoning model: from the ECG alone, it can identify several ECHO-derived structural and functional signals that are relevant for triage and follow-up imaging decisions. Echocardiography nevertheless remains necessary for definitive assessment of cardiac structure, ventricular function, and valve disease.

2.4 ECG-LLM answers clinically relevant questions across diverse ECG concepts

We next tested whether ECG-LLM understands a broad range of ECG concepts and can answer clinically relevant questions about raw ECG recordings. ECG-QA is a benchmark for question answering over ECG signals, built from ECG-expert-validated templates covering clinically relevant findings [37]. The benchmark includes questions about individual ECGs and comparisons between ECG recordings, with answer formats that test verification, choice between candidate findings, and open recovery of relevant findings. We evaluated ECG-LLM on the single-ECG subset from PTB-XL and MIMIC-IV-ECG, focusing on Single-Verify, Single-Choose, and Single-Query questions. These formats reflect increasing difficulty: Single-Verify asks whether a specific finding is present, Single-Choose requires discrimination between candidate findings, and Single-Query requires recovering the full set of relevant findings for a broader question. From a first-line clinical perspective, Single-Verify and Single-Choose are closest to the focused

questions a general practitioner or triage clinician may ask about a suspected ECG abnormality. Single-Query is still useful as a benchmark of broad, report-like recovery, but it is less central to our intended use case because it asks the model to enumerate a complete multi-label answer set rather than answer a targeted clinical question. The evaluated concepts include, but are not limited to, rhythm disorders, conduction disease, ischemic or infarction-related changes, hypertrophy, chamber enlargement, repolarization abnormalities, axis deviation, pacing, ectopy, interval abnormalities, and signal quality. Although ECG machines already provide intervals, axes, and rule-based diagnostic statements, these concepts remain important in ECG-QA because they test whether the model has learned the basic ECG descriptors required for broader ECG understanding. During evaluation, ECG-LLM received the question and candidate options, generated a free-text answer, and this answer was parsed back into the ECG-QA option space using regular-expression matching.

PTB-XL					MIMIC-IV-ECG				
Type	EM	μ F1	MF1	Maj	Type	EM	μ F1	MF1	Maj
Overall	57.71	65.60	67.86	39.28	Overall	43.75	59.36	64.10	24.83
Single-Verify	76.71	76.71	76.71	67.73	Single-Verify	81.06	81.06	81.06	66.45
Single-Choose	68.32	74.08	74.70	31.15	Single-Choose	68.37	73.49	74.57	30.91
Single-Query	38.26	58.35	57.75	23.19	Single-Query	19.69	53.70	53.39	5.61

Table 3: ECG-QA performance on PTB-XL and MIMIC-IV-ECG test sets. EM is exact-match accuracy on the set of selected options, μ F1 is micro-averaged F1 over all question-option decisions, MF1 is macro-averaged F1, and Maj is a majority baseline that always predicts the most frequent option set for each question type.

PTB-XL					
Type	Majority	Fusion [37]	MedViLL [38]	M3AE [39]	ECG-LLM
Single-Verify	67.73	72.10	73.90	74.60	76.71
Single-Choose	31.15	49.20	54.40	57.10	68.32
Single-Query	23.19	32.00	38.20	41.00	38.26

Table 4: Exact-match comparison on PTB-XL Single-Verify, Single-Choose, and Single-Query. Baseline values (Fusion, MedViLL [38], M3AE [39]) are taken from the ECG-QA benchmark and correspond to supervised multi-label classifiers over a fixed ECG-QA answer set. ECG-LLM instead generates free-text answers, which are parsed into label sets and evaluated with the same exact-match metric. Best score in each row is bolded.

ECG-LLM performed consistently above the per-question-type majority baseline on both datasets (Table 3). Overall, exact match improved from the majority baseline of 39 to 58 EM on PTB-XL and from 25 to 44 EM on MIMIC-IV-ECG. Unless otherwise stated, paired values report PTB-XL/MIMIC-IV-ECG. Option-level F1 was higher than exact match, indicating that many errors were partial rather than completely wrong interpretations (66/59 μ F1). Performance was strongest for targeted questions: Single-Verify reached 77/81 EM and Single-Choose reached 68/68 EM. Compared with published PTB-XL ECG-QA baselines, ECG-LLM was substantially stronger on the targeted ECG-QA tasks despite using a generative formulation rather than supervised classification over a fixed answer space (Table 4). It achieved the best reported result on Single-Verify (76.71 vs. 74.60 EM for M3AE [39]) and Single-Choose (68.32 vs. 57.10 EM for M3AE), giving the strongest overall performance on the clinically focused ECG-QA formats. On the harder Single-Query task, ECG-LLM ranked second (38.26 vs. 41.00 EM for M3AE) while still exceeding Fusion [37] and MedViLL [38]. The main remaining gap

was Single-Query exact match, where the model must recover a complete multi-label answer set. This metric is strict for free-text generation: missing one correct label or adding one extra label in the free-text answer makes the answer incorrect. Accordingly, Single-Query option-level F1 was substantially higher than exact match on both datasets. On PTB-XL, Single-Query reached 38 EM and 58 μ F1. On MIMIC-IV-ECG, it reached 20 EM and 54 μ F1, showing that many errors reflected incomplete answer sets rather than entirely incorrect ECG interpretations.

Stratifying the results by ECG-QA attribute type shows that the strong targeted-question performance is not driven by one narrow label group, but extends across many clinically interpretable ECG concepts (Supplementary Tables 10–11). In Single-Verify questions, exact match reached 68–90 EM across heart axis, infarction stage, extra systoles, SCP diagnostic statements, numeric features, and noise-related questions on PTB-XL, and 78–82 EM across all corresponding MIMIC-IV-ECG groups available in the benchmark (SCP diagnostic statements, heart axis, infarction stage, noise-related questions, and numeric features). In Single-Choose questions, ECG-LLM answered multiple concept groups: SCP-code diagnostic statements reached 69/68 EM (75/74 μ F1), heart-axis questions 75/50 EM, infarction-stage questions 56/78 EM, and numeric-feature questions 61/61 EM. Noise and extra-systole questions were lower on PTB-XL (43 and 50 EM).

For individual findings (Figures 4 and 5, Supplementary Table 12), we focused on clinically salient Single-Verify attributes because these most directly correspond to targeted yes/no questions about a suspected ECG finding. ECG-LLM was strong for major conduction disease, including complete left and right bundle branch block (100/93 and 97/98 EM, respectively), incomplete left bundle branch block (92/92 EM), left posterior fascicular block (98/85 EM), and first-degree AV block (67/83 EM). It also performed well for rhythm and pacing findings, including atrial fibrillation (87/90 EM), sinus tachycardia (98/93 EM), sinus bradycardia (83/95 EM), normal functioning artificial pacemaker (100 EM on PTB-XL), pacemaker rhythm (97 EM on MIMIC-IV-ECG), and ventricular-paced complexes or rhythm (95 EM on MIMIC-IV-ECG). Structural, infarction, and repolarization concepts were also recovered, including hypertrophy (88/83 EM), right atrial enlargement (91/85 EM), long QT interval (79/88 EM), T-wave abnormality (81/80 EM), Q waves (83/80 EM), and myocardial infarction as a broad diagnostic concept (78/85 EM). Additional high-performing clinically important labels included atrial flutter (94 EM on PTB-XL), atrial tachycardia (92 EM on MIMIC-IV-ECG), low QRS voltage (94 EM on MIMIC-IV-ECG), biventricular hypertrophy (92 EM on MIMIC-IV-ECG), and rapid ventricular-response or ventricular tachycardia-related labels such as non-sustained ventricular tachycardia, wide QRS tachycardia, and uncontrolled ventricular response (93–98 EM on MIMIC-IV-ECG). Axis and signal-quality questions were also clinically useful, with axis-deviation findings reaching 67–93 EM and noise-related Single-Verify performance reaching 68/79 EM at the group level. The main limitations were more granular tasks. Localized infarction or ischemia labels, which require assigning abnormalities to specific lead territories such as anterior, inferior, lateral, or septal regions, generally ranged from 58–86 EM. Left atrial overload/enlargement remained difficult (55/44 EM). Interval-range questions varied (62–98 EM on PTB-XL and 67–100 EM on MIMIC-IV-ECG). Broad Single-Query questions remained harder (38/20 EM, despite 58/54 μ F1) because exact match penalizes incomplete multi-label answers. Overall, this ECG-QA evaluation shows that ECG-LLM has learned a broad ECG vocabulary and can answer focused, clinically relevant questions from the raw signal, while complete report-like recovery, fine-grained localization, subtle atrial findings, and exact interval reasoning remain important limitations.

2.5 ECG-LLM generates free-text ECG diagnostic summaries

We next evaluated whether ECG-LLM can generate free-text ECG diagnostic summaries from raw ECGs. ECG machines already provide automated measurements and diagnostic statements, and many prior ECG-language and ECG-image models use report generation or diag-

ECG-LLM performance on ECG-QA Benchmark on PTB-XL

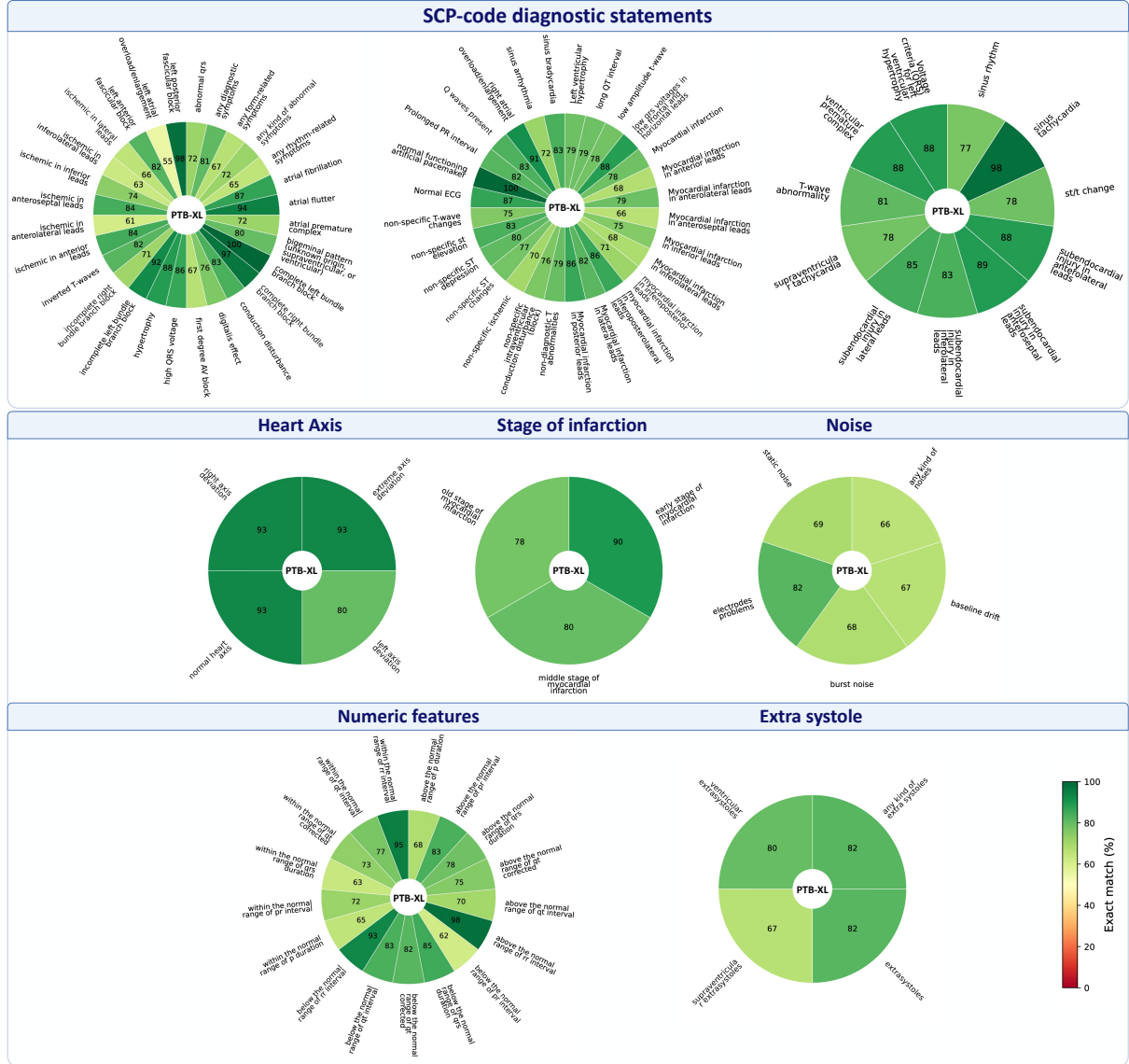


Figure 4: Single-Verify ECG-QA exact-match (EM) performance across ECG concepts on PTB-XL. Results include SCP-code diagnostic statements, heart axis, stage of infarction, noise, numeric features, and extra systoles. Sector values indicate exact-match performance for each attribute. The full table is provided in the Supplementary Table 12.

ECG-LLM performance on ECG-QA Benchmark on MIMIC-IV-ECG

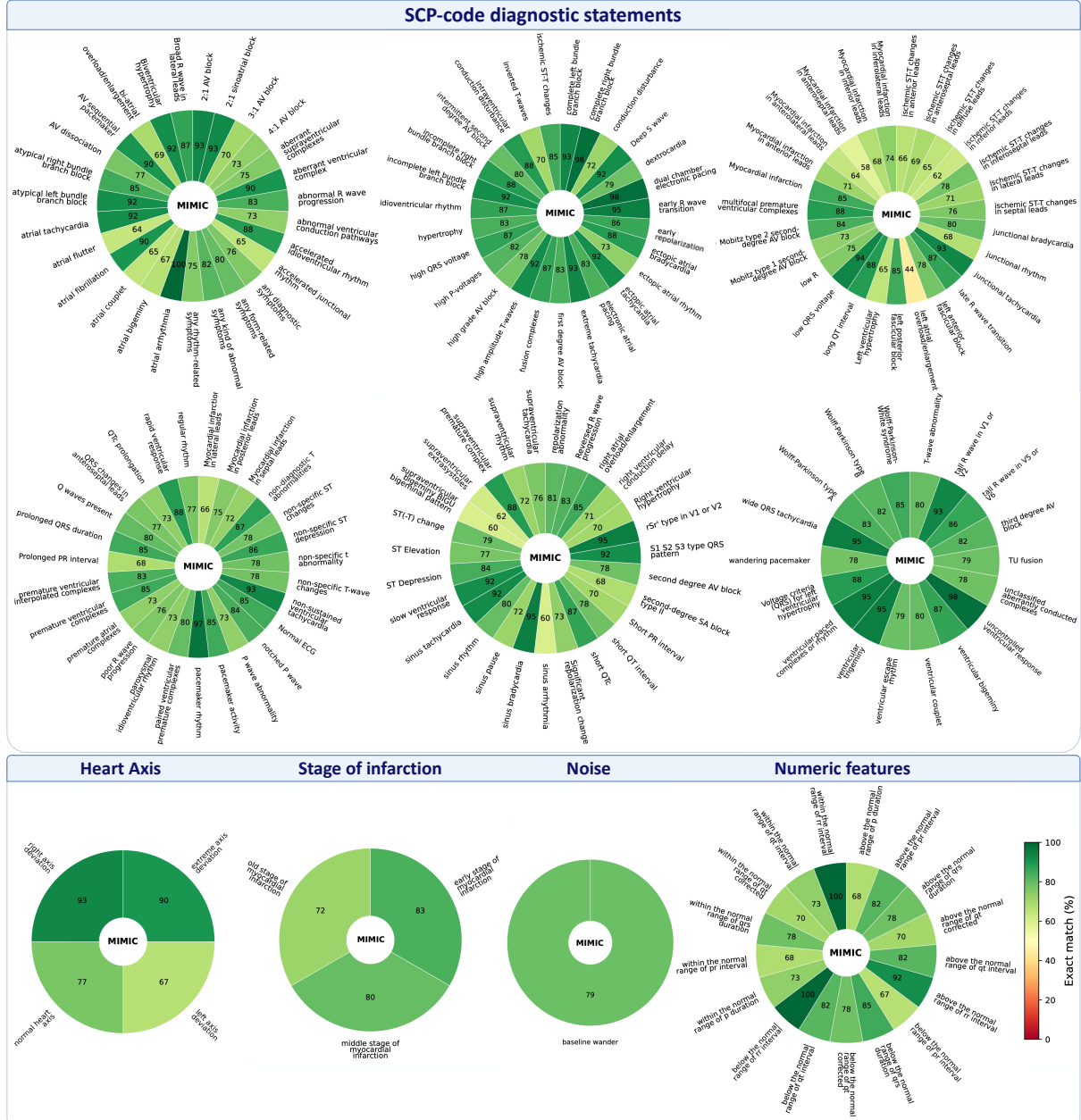


Figure 5: Single-Verify ECG-QA exact-match (EM) performance across ECG concepts on MIMIC-IV-ECG. Results include SCP-code diagnostic statements, heart axis, stage of infarction, noise, and numeric features. Sector values indicate exact-match performance for each attribute. The full table is provided in the Supplementary Table 12.

nostic interpretation as a primary endpoint. In our study, this task serves a different purpose: it tests whether the ECG-conditioned language model has learned the mapping between raw ECG morphology and conventional diagnostic language. Diagnostic reports provide dense weak supervision because they summarize rhythm, conduction, repolarisation, ischemia/infarction, pacing, and other ECG abnormalities in compact clinical language. Solving this task, therefore, demonstrates that the model can associate ECG signal patterns with clinically named ECG concepts, strengthening the ECG-language grounding needed for broader question answering over ECG-derived and imaging-related cardiac phenotypes. For each ECG and each evaluated model, one diagnostic-report question was randomly selected from a fixed prompt set (Supplementary Table 17), and the model was asked to produce a concise free-text ECG diagnosis. The generated report was compared with the reference diagnosis by an LLM judge scoring correctness, completeness, and specificity.

We sampled 1,000 non-trivial ECG studies from unique patients in each of UK Biobank [10], PTB-XL [12], and MIMIC-IV-ECG [11]. Non-trivial records were defined using dataset-specific normal-report rules. In the UK Biobank, an ECG was considered non-trivial if the device diagnosis was not exactly the pair "Normal sinus rhythm" and "Normal ECG". In MIMIC-IV-ECG, an ECG was considered non-trivial if the device diagnosis was not exactly "Sinus rhythm; Normal ECG". In PTB-XL, an ECG was considered non-trivial if its SCP-code set was not exactly {NORM, SR}. For each selected ECG, ECG-LLM was asked to produce a concise free-text ECG diagnosis. Ground truth was obtained from dataset-specific diagnostic annotations: device diagnoses for UK Biobank and MIMIC-IV-ECG, and report-derived diagnostic annotations for PTB-XL (the original German/Swedish report and its English translation).

We compared ECG-LLM with three publicly available baseline models: MedGemma-4B [40, 41, 42, 43], ECG-Chat [33], and PULSE-7B [34, 44]. ECG-Chat is the closest raw-signal ECG-language baseline. It’s published system combines an ECG encoder and LLM with additional components for report generation, including diagnosis-driven prompting, GraphRAG-based retrieval [45], DSPy-based prompt tuning [46], and a structured LaTeX report-generation pipeline [33]. Here, we evaluate ECG-Chat without GraphRAG retrieval or downstream report-generation modules to focus on ECG-conditioned generation rather than a retrieval-augmented reporting system, although results with GraphRAG may be higher because it uses an additional supervised ECG interpretation method. PULSE-7B is a strong ECG-image baseline trained for ECG image interpretation using ECGInstruct, a large-scale ECG image-text instruction dataset [34, 44, 47]. For this comparison, raw ECG signals were rendered as ECG images using ECG-Image-Kit [48, 49]. UK Biobank is not among the reported training sources of PULSE-7B, so its UK Biobank evaluation should be interpreted as out-of-domain. MedGemma-4B was included to show the current performance of a general medical vision-language model on the same rendered ECG images. Although MedGemma-4B has been trained with medical image data that includes ECG images, it is not an ECG-specialised model like PULSE-7B. Thus, the baselines cover both raw-signal and image-based ECG-language settings, while all models are assessed with the same diagnostic-generation prompt, reference reports, and LLM-judge protocol.

Standard NLP string-overlap metrics provide an incomplete evaluation of free-text ECG reports. The same clinical finding can be expressed with different wording, while similar wording can correspond to different clinical meanings. Moreover, the severity of an error depends on the finding, so missing a major rhythm, conduction, or ischemic abnormality is not equivalent to omitting a minor or borderline statement. We therefore used a locally deployed GPT-OSS-20B model [50, 51, 52] as a reasoning-capable LLM judge, following the broader use of LLM-as-judge evaluation for open-ended model outputs [53] and recent work applying LLM judges to medical report generation [54]. The judge acted as an automated clinical adjudicator: it compared the generated report with the reference diagnosis, assigned structured scores for correctness, completeness, and specificity, and provided notes as its reasoning. Correctness measures whether

LLM-as-Judge Evaluation of free-text ECG report generation

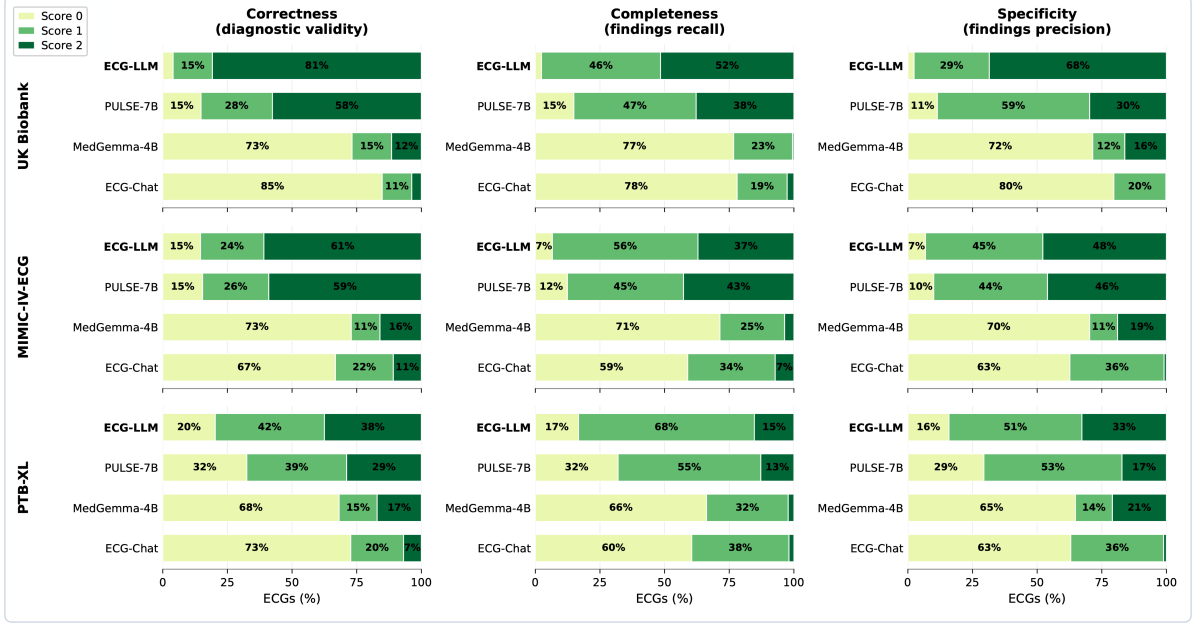


Figure 6: LLM-judge evaluation of free-text ECG diagnostic summaries across UK Biobank, MIMIC-IV-ECG, and PTB-XL. A locally deployed GPT-OSS-20B judge compares each generated diagnosis with the ground-truth diagnostic report and scores correctness, completeness, and specificity on a three-level scale. Score 0 denotes poor performance, score 1 partial performance, and score 2 good performance. Bars show the percentage of ECGs assigned to each score level.

the main ECG interpretation is clinically accurate overall, completeness measures whether important ground-truth findings are recovered, and specificity measures whether the generated findings are supported by or compatible with the reference report. Each score ranges from 0 to 2, where 0 indicates poor performance, 1 indicates partial performance, and 2 indicates good performance. The judge’s prompt allowed clinically equivalent wording while penalising major rhythm, conduction, ischemia/infarction, pacing, and unsupported-diagnosis errors. The full judge prompt is provided in Supplementary Section 15.

Figure 6 shows that ECG-LLM and PULSE-7B are the strongest models in this evaluation. On UK Biobank, ECG-LLM achieves the clearest advantage, with a substantially larger fraction of score-2 reports for correctness and specificity. PULSE-7B nevertheless remains a strong baseline, especially because UK Biobank is out-of-domain for this model and because it receives the ECG only as a rendered image. On MIMIC-IV-ECG, ECG-LLM and PULSE-7B are broadly comparable: both models obtain similar score distributions for correctness and specificity, while PULSE-7B shows slightly higher completeness. On PTB-XL, ECG-LLM shows higher fractions of score-2 correctness and specificity than PULSE-7B, although both models frequently receive partial completeness scores, indicating that open-ended reports often recover some but not all reference findings. For PTB-XL, the LLM judge used the original PTB-XL report as the primary reference. The English translation was provided only to help interpret the non-English report text and was not treated as an additional list of required findings, as specified in the judge prompt. MedGemma-4B and ECG-Chat perform less consistently under this protocol. For MedGemma-4B, this likely reflects the difficulty of ECG interpretation for a general medical vision-language model, even when ECG images are included in its broader medical-image training data. For ECG-Chat, the result should be interpreted as a model-level raw-signal comparison rather than as an evaluation of its full published ECG-Chat system. Its published report-generation pipeline explicitly supplements the ECG-conditioned model

with diagnosis-driven prompting, retrieval-augmented knowledge with GraphRAG, prompt optimization with DSPY, and structured report formatting. Therefore, the lower performance observed here is consistent with the smaller scale and narrower coverage of the underlying ECG-language instruction data, and with the possibility that part of ECG-Chat’s long-form reporting capability depends on its retrieval and report-generation modules rather than on the bare ECG-conditioned generator alone.

Overall, these results show that ECG-LLM learns a useful mapping from raw ECG morphology to conventional ECG diagnostic language. This report-generation task is not the main clinical aim of ECG-LLM: unlike prior ECG-language systems that focus primarily on reproducing diagnostic reports, our central goal is flexible question answering over ECG-derived and imaging-related cardiac phenotypes. Nevertheless, strong performance on diagnostic summaries is an important intermediate learning step, because it confirms that the ECG-conditioned language model grounds signal patterns in clinically named ECG concepts. ECG-LLM remains competitive with the strongest ECG-image baseline and outperforms the general medical image-text baseline and the raw-signal ECG-language baseline in most settings, while the partial completeness scores observed across models highlight the continuing difficulty of open-ended report generation, especially for secondary findings and dataset-specific report annotations.

2.6 Qualitative evaluation of ECG-LLM

We assessed whether ECG-LLM produces ECG-conditioned outputs while preserving general language competence. When no ECG waveform was provided, the model declined to make a diagnosis and requested the missing trace. When the ECG input was replaced by random noise, the output lost clinical coherence, suggesting that meaningful predictions depend on structured ECG information rather than language priors alone. For non-cardiac general-knowledge prompts without ECG tokens, the model behaved similarly to the underlying instruction-tuned LLM, indicating no apparent catastrophic forgetting after ECG-specific fine-tuning.

We present qualitative examples from UK Biobank [10], EchoNext [13], MIMIC-IV-ECG [11], and PTB-XL [12] (Figures 7–9). The UK Biobank and EchoNext examples include both ECG-derived findings, such as rhythm, intervals, axis, and conduction abnormalities, and ECG-predicted imaging-derived phenotypes, including ventricular function, chamber volumes, valvular disease, and pulmonary pressure-related findings. The MIMIC-IV-ECG and PTB-XL examples focus on direct ECG interpretation, including rhythm, interval measurements, conduction abnormalities, ischemic changes, and overall diagnostic impressions.

3 Discussion

This work introduces ECG-LLM, an ECG-conditioned language model for flexible cardiac question answering from a 12-lead ECG. Our central contribution is to use language as a common supervision space for heterogeneous cardiovascular observations. ECG measurements and reports, CMR-derived phenotypes, echocardiographic phenotypes, and diagnostic labels can all be expressed as clinically grounded question–answer pairs and used to train a single ECG-conditioned model. The work is motivated by broader shifts in ECG AI: from conventional ECG-derived interpretation toward language-based ECG assistants, and from ECG-only phenotype prediction toward cross-modal learning of cardiac phenotypes that are not directly measured on the ECG [33, 34, 35, 26, 27]. Our emphasis is to combine these directions. ECG-LLM uses question–answer supervision to connect raw ECG morphology with both conventional ECG concepts and imaging-derived structural and functional phenotypes, including ventricular function, chamber size, wall thickness, and valve- or pressure-related findings. In a front-line setting, this creates an ECG-only interface through which clinicians can ask focused questions about the ECG and broader cardiovascular state to support triage, follow-up imaging, or referral prioritization.

Qualitative Examples of ECG-LLM Responses on UK Biobank ECGs

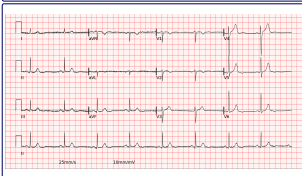
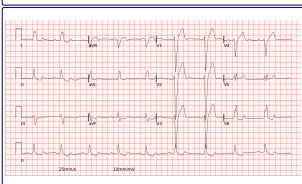
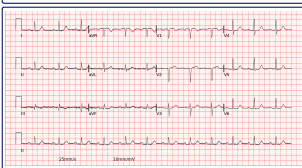
	Q1: What is the diagnostic summary of this ECG? A1: The ECG demonstrates a sinus rhythm with a left bundle branch block pattern. The QRS duration is prolonged, and the QTc interval is also extended. The T axis is rightward, while the R axis remains normal. Heart rate and RR interval are within normal limits, and the PQ interval is normal. GT: QTc = 475 [ms]; male Q2: Review the ECG and describe any abnormalities in LV structure and function. A2: The ECG shows a dilated left ventricular end systolic volume, a dilated left ventricular end diastolic volume, a mildly reduced left ventricular ejection fraction, and a borderline low cardiac index. GT: cardiac index = 2.9 [litres/min/m²]; male
	Q1: What abnormalities are present on this ECG? A1: The ECG demonstrates marked sinus bradycardia with a markedly prolonged RR interval and a prolonged PP interval, indicating a slow heart rate. The left ventricular wall thickness appears thinned, and the cardiac index is borderline low, pointing to reduced cardiac output. Overall, the tracing is abnormal. GT: cardiac index = 2.8 [litres/min/m²]; female Q2: Review this ECG and describe ventricular and atrial function and structure. A2: The ECG demonstrates a normal left ventricular ejection fraction, normal left ventricular end-diastolic and end-systolic volumes, normal left ventricular stroke volume, and normal left ventricular myocardial mass. The right ventricular ejection fraction, end-diastolic and end-systolic volumes, and stroke volume are all within normal limits. Left atrial and right atrial ejection fractions, stroke volume, and maximum and minimum volumes are normal. GT: LA min volume (indexed): 24 [mL/m²]; female; BSA = 1.47 [m²]
	Q1: Give a summary of ventricular function and structure. A1: The ECG reveals a dilated RV end systolic volume, indicating an abnormality in the right ventricular chamber. The ventricles exhibit a reduced right ventricular ejection fraction and a reduced left ventricular ejection fraction, with a dilated left ventricular end systolic volume. Furthermore, the RV stroke volume is increased, which may be a compensatory mechanism in response to the reduced ejection fraction. GT: RV EF = 44%; RV stroke volume = 128 [mL]; male Q2: Describe the heart rhythm, cardiac axis, intervals, and any conduction abnormalities on this ECG. A2: The ECG shows sinus bradycardia with a left bundle branch block, prolonged RR interval and a prolonged PP interval, a leftward T axis, a leftward P axis, and a prolonged QRS duration. Overall, the ECG is abnormal, suggesting potential cardiac conduction abnormalities. GT: P axis: 28 [degrees]; male
	Q1: Provide an ECG report for this tracing. A1: The ECG demonstrates sinus rhythm with a first-degree atrioventricular block, evidenced by a prolonged PQ interval. The P wave duration is also prolonged, while the QRS complex, QTc interval, and all electrical axes remain within normal limits. Heart rate and RR interval are normal. Additionally, the ECG suggests the possibility of a septal infarct, or a heart attack affecting the septum, although the age of the infarct cannot be determined. GT: P duration = 92 [ms]; Missed – "Possible Left atrial enlargement"; female

Figure 7: Qualitative ECG-LLM examples from UK Biobank for ECG-derived findings and imaging-derived cardiac phenotypes. Representative cases show the ECG input, the user question, the model answer, and the corresponding ground-truth label or measurement. Green highlights denote exact agreement with the ground truth, yellow highlights denote clinically close or partially correct predictions, and red highlights denote incorrect predictions. For non-exact predictions, the relevant ground-truth value or label is shown below the model answer.

This design has several practical advantages. With the same model, a single ECG representation can support open-ended ECG questions, imaging-related phenotype queries, benchmark-style QA, and report-style summaries without changing the architecture for each task. This is clinically attractive because the relevant question often depends on symptoms, context, and uncertainty, and may not be known when a model is designed. The language-supervision framework also supports scalability. Many cohorts may not be able to share raw CMR or echocardiography images, but structured reports, measurements, or extracted phenotypes can often be represented as text; once expressed in this form, they can enter the same QA-generation pipeline. Thus, ECG-LLM can be trained from ECG waveforms paired with text-derived supervision from multiple modalities, without requiring raw imaging data during instruction tuning. By including both ECG-derived and imaging-derived phenotypes in QA generation, the model is encouraged to learn associations between ECG morphology, conventional ECG interpretation, and broader cardiovascular status rather than only ECG-machine-style outputs. The ECG-QA benchmark and report-generation experiments therefore serve mainly as grounding checks, showing that the model has learned ECG concepts and a mapping between ECG signal patterns and conventional ECG language, while the broader aim is flexible second-opinion support for triage and referral prioritization.

Several limitations remain. First, ECG-LLM is still an LLM-based system, so hallucination and overconfident generation remain important risks. These risks can arise during free-text answering, structured extraction, and LLM-based evaluation, and they are especially important in clinical settings where plausible wording can mask incorrect content. Prior work has shown that medical LLMs can encode useful clinical knowledge while still producing factual errors and potentially harmful recommendations, with performance sensitive to prompting and input

Qualitative Examples of ECG-LLM Responses on EchoNext ECGs

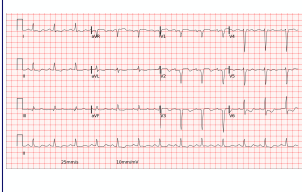
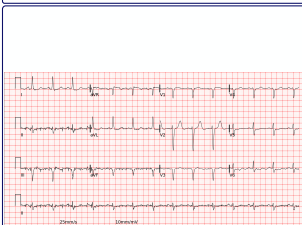
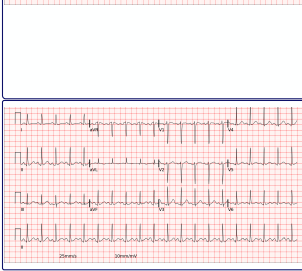
	Q1: Does this ECG suggest reduced left ventricular systolic function? A1: Yes, the ECG indicates severely reduced left ventricular systolic function, with an ejection fraction of 45% or lower. Q2: Does this ECG suggest moderate-to-severe right ventricular systolic dysfunction? A2: Yes, the ECG indicates that right ventricular systolic dysfunction is moderate or greater, and right ventricular systolic function is severely reduced. Q3: Does this ECG show moderate left ventricular hypertrophy with the maximum of the interventricular septum (IVS) or posterior wall (LVPW) thickness is greater than or equal to 1.3 cm? Think step by step and then give your final answer. A3: To determine if this ECG shows moderate left ventricular hypertrophy, we need to examine the thickness of the interventricular septum (IVS) and the left ventricular posterior wall (LVPW). The ECG indicates that the IVS thickness is mildly increased, which is a sign of left ventricular hypertrophy. However, the LVPW thickness is not mentioned as being greater than or equal to 1.3 cm. In fact, the LVPW thickness is reported as normal. Therefore, based on the available information, we cannot conclude that the ECG shows moderate left ventricular hypertrophy, as the LVPW thickness does not meet the specified criterion.
	Q1: Does this ECG suggest a moderate or large pericardial effusion? A1: The ECG indicates that there is no moderate or large pericardial effusion. Q2: Does this ECG suggest moderate or severe mitral regurgitation? A2: The ECG indicates that mitral regurgitation is moderate or greater, and it is specifically classified as severe. This reflects a significant regurgitant flow across the mitral valve. Q3: Does this ECG show moderate or severe aortic stenosis? A3: The ECG indicates that aortic stenosis is moderate or greater, and it is classified as severe. This is the most prominent valvular abnormality noted on the ECG. GT: Aortic stenosis = "moderate", female Q4: Does this ECG suggest moderately or severely reduced left and/or right ventricular systolic function? Think step by step and give your final answer. A4: Let's break down the information step by step: 1. The ECG indicates that the left ventricular systolic function is severely reduced, as it meets the criteria for left ventricular ejection fraction of 45% or lower. 2. The right ventricular systolic function is also severely reduced, as it meets the criteria for right ventricular systolic dysfunction of moderate or greater severity. 3. The presence of severely reduced left and right ventricular systolic function suggests that both ventricles are experiencing significant impairment in their pumping abilities. 4. The combination of these findings indicates that the ECG suggests moderately or severely reduced left and right ventricular systolic function. Final answer: The ECG suggests moderately or severely reduced left and right ventricular systolic function. GT: Left ventricular systolic function = 37.5%, "moderately reduced"; Right ventricular systolic function = "moderately reduced"; female
	Q1: Does this ECG suggest pulmonary hypertension, as would be expected with a pulmonary artery systolic pressure of at least 45 mmHg? Think step by step and give your final answer. A1: Let's break down the information step by step. First, we need to identify any signs of pulmonary hypertension on the ECG. The ECG does not show any right ventricular hypertrophy or strain pattern, which are common indicators of pulmonary hypertension. However, the ECG does indicate that the pulmonary artery systolic pressure is at least 45 mmHg, which is a key indicator of pulmonary hypertension. Given the absence of right ventricular hypertrophy or strain pattern, it is likely that the pulmonary hypertension is not severe or is compensated. Therefore, the ECG suggests pulmonary hypertension, but it is not severe. Q2: Does this ECG suggest moderately or severely reduced left ventricular systolic function? A2: Left ventricular systolic function is normal, and the ejection fraction is not 45% or lower, so there is no evidence of moderate or severe systolic dysfunction.

Figure 8: Qualitative ECG-LLM examples from EchoNext for ECG-derived findings and imaging-derived cardiac phenotypes. Representative cases show the ECG input, the user question, the model answer, and, where applicable, the corresponding ground-truth diagnosis or measurement. Green highlights denote exact agreement with the ground truth, yellow highlights denote clinically close or partially correct predictions, and red highlights denote incorrect predictions. For non-exact predictions, the relevant ground-truth value or label is shown below the model answer.

formulation [55]. ECG-LLM should therefore be viewed as clinician-facing decision support, not as an autonomous diagnostic system or a patient-facing replacement for expert interpretation. Clinical human evaluation also remains necessary before deployment. We partly addressed the absence of large-scale expert review by using an LLM judge for scalable report-generation assessment, but this should be viewed as a screening and comparison tool rather than a substitute for cardiologist evaluation, especially for safety-critical errors, subtle findings, and clinical actionability.

Second, some evaluated phenotypes are continuous measurements that were converted into discrete categories. This approach facilitates evaluation and is better suited to language-model supervision, as textual categories are generally more natural targets for generative models than exact numerical regression. However, discretization introduces boundary effects, whereby values close to a threshold may be clinically borderline but still counted as categorical errors. The free-text formulation partly alleviates this issue by allowing the model to describe a finding as borderline or uncertain, although future work should more directly incorporate continuous uncertainty, calibrated numerical estimates, and clinically meaningful threshold handling. Dataset imbalance represents an additional limitation, particularly in UK Biobank [10], where many cardiac phenotypes are dominated by normal values. This imbalance can reduce the representation of abnormal findings during training and evaluation. We partly mitigate this imbalance by incorporating datasets with a higher prevalence of abnormal findings and by enriching the UK Biobank training data during QA generation. For studies containing at least one abnormal field, we retained the available abnormal findings and generated additional QA pairs to increase the model's exposure to non-normal phenotypes.

Qualitative Examples of ECG-LLM Responses on MIMIC-IV-ECG and PTB-XL ECGs

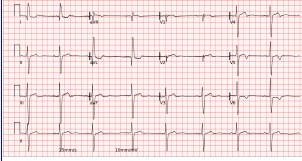
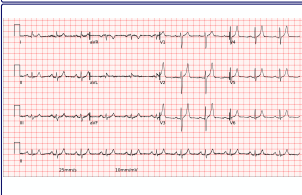
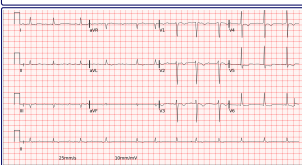
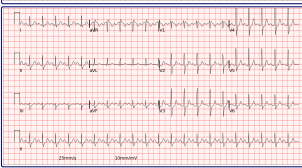
	Q1: Review this ECG signal and produce a detailed report based on your diagnostic observations. A1: The electrocardiogram (ECG) reveals a sinus bradycardia with sinus arrhythmia, with a left axis deviation suggesting an abnormal direction of the heart's electrical activity. Additionally, the ECG shows signs of a right bundle branch block (RBBB) with a left anterior fascicular block, which may indicate a conduction abnormality in the heart's electrical system. Furthermore, the ECG suggests possible evidence of a lateral infarct, or myocardial infarction, of undetermined age, as well as possible left ventricular hypertrophy, which may be contributing to the inferior and lateral ST-T changes observed. GT: ECG diagnosis = Sinus bradycardia with sinus arrhythmia; Left axis deviation; Iv conduction defect; Lateral infarct - age undetermined; Left ventricular hypertrophy; Abnormal ECG; QRS = 136 [ms]; R axis = -57°; female Q2: Is the QRS prolonged on this ECG or not? A2: Yes, the QRS duration is prolonged, which is a key feature of the intraventricular conduction defect.
	Q1: Does this ECG show any evidence of arrhythmia, conduction disease, or acute myocardial ischemia? A1: Let me think through this step by step: Step 1: I first examine the rhythm and rate. The ECG diagnosis states normal sinus rhythm, and the heart rate is normal. Step 2: I then look at conduction intervals: the PQ interval, QRS duration, and QTc interval are all reported as normal. Step 3: I also check for any abnormal axis or interval prolongation that could suggest conduction disease. The P axis, R axis, and T axis are all normal, and the QRS duration is normal. Step 4: Finally, I review the ECG for any ischemic changes such as ST-segment elevation or depression, but the ECG is described as normal. Therefore, there is no evidence of arrhythmia, conduction disease, or acute ischemia. So the answer is: No, the ECG shows no evidence of arrhythmia, conduction disease, or acute myocardial ischemia. Q2: Review the ECG signal and produce a detailed report on your diagnostic observations, ending with the final diagnosis. A2: The electrocardiogram (ECG) reveals a normal sinus rhythm, indicating that the heart is beating at a regular rate and rhythm, with a normal P-wave, PR interval, QRS complex, and QT interval. The ECG tracing does not show any signs of arrhythmia, conduction abnormalities, or ischemic changes. The overall impression is that the ECG is normal, suggesting that the heart's electrical activity is functioning within normal limits. Final Diagnosis: Normal ECG.
	Q1: In reference to this ECG, can you determine which diagnostic symptom is being shown, ischemic in anteroseptal leads or myocardial infarction in inferior leads, even if they are not entirely clear? A1: The ECG suggests the presence of ischemic changes in the anteroseptal leads, which may be indicative of myocardial ischemia in this region. Additionally, there is a possibility of myocardial infarction in the inferior leads, although this is not entirely clear. The combination of these findings may indicate a more complex cardiac condition, warranting further evaluation and diagnosis. Q2: Does this ECG reveal any evidence of atrial fibrillation? A2: Yes, this ECG reveals evidence of atrial fibrillation. The presence of atrial fibrillation is indicated by the irregular rhythm, which is characterized by the absence of P waves and a lack of a consistent atrial activity.
	Q1: Which ECG measurements have values within the normal range? Possible options: P duration, PR interval, QRS, QTc, QT interval or RR interval. A1: qrs duration, qt corrected, p duration Q2: What is heart rhythm observed on this ECG? A2: The heart rhythm observed on this ECG is sinus tachycardia, which is a type of abnormal heart rhythm characterized by a rapid heart rate. Q3: Can you see incomplete right bundle branch block on this ECG? A3: Yes, the ECG shows an incomplete right bundle branch block. This is evident from the characteristic slurred S wave in lead V1, although the overall QRS duration remains within normal limits.

Figure 9: Qualitative ECG-LLM examples from MIMIC-IV-ECG and PTB-XL for direct ECG interpretation. Representative cases show the ECG input, the user question, the model answer, and, where applicable, the corresponding ground-truth diagnosis or measurement. The examples cover diagnostic report generation and questions related to rhythm, interval measurements, conduction abnormalities, and ischemic findings. Green highlights denote exact agreement with the ground truth, yellow highlights denote clinically close or partially correct predictions, and red highlights denote incorrect predictions. For non-exact predictions, the relevant ground-truth value or label is shown below the model answer.

Third, ECG-LLM infers imaging-derived cardiovascular phenotypes from patterns represented in the ECG rather than directly measuring cardiac structure or function. Previous studies have demonstrated that ECGs contain predictive information about ECHO- and CMR-derived structural and functional phenotypes, supporting the feasibility of learning cross-modal cardiovascular information from ECG signals [13, 26]. Nevertheless, the strength of this association varies across phenotypes, and some imaging-derived findings may be more readily inferred from ECG morphology than others.

Architecturally, ECG-LLM currently passes all ECG encoder tokens to the language model through the projector, allowing the model to retain the full encoded ECG representation. Future work could explore more compact conditioning strategies, such as learnable pooling or attention-based token selection, to summarize the ECG sequence into a smaller set of task-relevant latent tokens before passing it to the LLM. Such approaches may improve computational efficiency and reduce context length while maintaining the diagnostically relevant information contained in the ECG.

Overall, ECG-LLM shows that knowledge derived from multiple cardiac modalities can be used as supervision for an ECG-conditioned question-answering model. By representing cardiovascular knowledge from imaging and non-imaging tests as language supervision, the approach turns the widely available 12-lead ECG into a more flexible interface for front-line cardiovascular reasoning when imaging or specialist review is delayed or unavailable. This expands the role of ECG-language modelling from generating reports toward answering clinically meaningful questions about a patient’s broader cardiovascular status.

4 Methods

4.1 Overview

ECG-LLM follows the common vision–language model design used in systems such as MiniGPT-4 and LLaVA [56, 57, 58]: a pretrained signal encoder extracts modality-specific tokens, a lightweight projector maps these tokens into the language-model embedding space, and an autoregressive LLM generates the text output. In our case, the visual encoder is replaced by a self-supervised ECG encoder, and the resulting model is trained to answer cardiac questions conditioned on the ECG. More specifically, given a 12-lead ECG $x \in \mathbb{R}^{C \times T}$ (with $C = 12$ leads and T time samples) and a question q , the ECG encoder produces latent representations $h = f_{\text{ecg}}(x)$, which are mapped into the LLM embedding space by a projector $z = g_{\text{proj}}(h)$. The LLM f_{llm} then generates an answer token sequence $y = (y_1, \dots, y_N)$ autoregressively conditioned on (z, q) . We train the model using the standard causal language modelling (next-token prediction) loss on the ground-truth answer tokens:

$$\mathcal{L} = - \sum_{n=1}^N \log p(y_n | y_{<n}, q, z), \quad z = g_{\text{proj}}(f_{\text{ecg}}(x)). \quad (1)$$

Our aim is to study how much clinically relevant information can be transferred into an ECG-conditioned language model through multimodal language supervision, while keeping inference strictly ECG-only. We therefore avoid task-specific ECG heads, retrieval pipelines, and prompt augmentation with precomputed ECG measurements at test time. Our model receives only the raw 12-lead ECG and the user’s question. This design allows us to isolate what can be inferred directly from the ECG representation under a unified free-text question-answering setup.

4.2 Architecture

ECG encoder. We use a self-supervised ECG encoder f_{ecg} pretrained with a masked autoencoder (MAE) objective [59]. The encoder follows a ViT-MAE design adapted to 12-lead ECGs. A patch embedding layer tokenizes each ECG into non-overlapping temporal patches, learnable positional embeddings are added, and a random subset of tokens is masked during pretraining. A Transformer [60] encoder processes the visible tokens, and a lightweight Transformer decoder is trained to reconstruct the masked patches. Let $\mathcal{M}$ denote the set of masked token positions, x_i the original ECG patch at masked position i , and $\hat{x}_i$ the corresponding reconstructed patch. The MAE reconstruction loss is:

$$\mathcal{L}_{\text{MAE}} = \frac{1}{|\mathcal{M}|} \sum_{i \in \mathcal{M}} \|\hat{x}_i - x_i\|_2^2. \quad (2)$$

Here, $|\mathcal{M}|$ is the number of masked tokens, and the loss is computed only over masked patches. After pretraining, we discard the decoder and retain only the encoder for ECG-LLM training.

We use MAE pretraining rather than ECG–text or ECG–imaging contrastive pretraining to keep the ECG encoder general-purpose and independent of paired reports or imaging data. This objective is well suited to ECGs because 12-lead traces contain local temporal structure and inter-lead redundancy, allowing masked segments to be reconstructed from surrounding waveform information.

Projector. The projector g_{proj} maps ECG encoder tokens into the LLM input embedding space. We write it as:

$$z = g_{\text{proj}}(h) = g_{\text{map}}(g_{\text{agg}}(h)), \quad (3)$$

where g_{agg} controls token aggregation and g_{map} performs the embedding-space mapping. In this work, we use identity aggregation, so all ECG tokens are preserved. The mapping module is a token-wise two-layer MLP with an activation function between the two linear layers.

Language model. The language component of ECG-LLM is an autoregressive LLM, Llama-3.1-8B-Instruct [61, 62]. The LLM receives the projected ECG tokens together with the tokenized question and predicts the answer tokens using the causal language-modelling objective (Eq. 1). Each example is formatted as a chat conversation: a fixed system prompt defines the model as a cardiology assistant, the user message begins with "Here is the ECG:" followed by the ECG tokens and the clinical question, and the assistant message contains the target answer. The full chat template is provided in the Supplementary Section 11. During ECG-LLM training, the original LLM parameters are kept frozen, and only LoRA adapters [63] are trained in selected self-attention projections. This reduces the risk of catastrophic forgetting and helps preserve the model’s broad general knowledge while still allowing adaptation to ECG-conditioned question answering. All implementation details are provided in Section 4.5 and Supplementary Section 13.

4.3 ECG cohorts and preprocessing

We use four 12-lead ECG cohorts to train and evaluate ECG-LLM: UK Biobank [10], MIMIC-IV-ECG [11], PTB-XL [12], and EchoNext [13]. These cohorts were chosen because they provide complementary views of cardiac status. MIMIC-IV-ECG and PTB-XL support conventional ECG interpretation through measurements, diagnostic statements, reports, and structured diagnostic labels. UK Biobank links ECGs to CMR-derived cardiac structure and function. EchoNext links ECGs to echocardiographic measurements and structural heart disease labels. Together, they allow the model to learn from both ECG-visible findings and imaging-related cardiac phenotypes. In total, the study includes 918,819 ECG recordings from 251,492 unique patients. The training split contains 679,112 ECGs from 186,409 patients. The validation split contains 68,151 ECGs from 18,615 patients. The test split contains 171,556 ECGs from 46,468 patients. Dataset split statistics are provided in Supplementary Table 5. For MIMIC-IV-ECG and PTB-XL, we adopt the official ECG-QA [37] splits to preserve benchmark comparability. For EchoNext, we use the official dataset splits. For UK Biobank, we construct a subject-wise split with 85% of patients assigned to the training set and 15% to the test set, which avoids patient overlap between training and testing.

For all datasets, we harmonise ECG waveforms to a shared 12-lead representation for the common ECG encoder. UK Biobank, MIMIC-IV-ECG, and PTB-XL waveforms are retained in physical millivolt units, while EchoNext waveforms are first denormalised using the released per-lead normalisation parameters and converted back to millivolts. We then reorder all ECGs to the PTB-XL lead order [I, II, III, aVL, aVR, aVF, V1–V6]. EchoNext ECGs are additionally resampled from 250 Hz to 500 Hz to match the common temporal grid. When baseline correction has not already been applied by the acquisition system, slow non-physiological baseline drift is removed lead-wise using a level-10 Daubechies-2 discrete wavelet transform by discarding the approximation coefficients before signal reconstruction, thereby reducing baseline wander while preserving diagnostically relevant ECG morphology [64].

4.4 Clinical question-answer supervision from cardiac records

To cover diverse cardiac concepts and support targeted questions that arise in general practice or front-line care, we use question-answer (QA) supervision. This format expresses dense cardiac information in natural clinical language and matches how clinicians ask patient-specific questions. QA pairs are also a standard format for instruction tuning multimodal assistants [57, 30]. In our setting, this strategy is useful because large-scale ECG datasets rarely provide diverse text instructions for question-driven cardiac reasoning.

We therefore construct ECG-specific instruction-tuning data from structured cardiac records (Figure 10). We refer to these generated examples as structured-record QA, because the answer source is a per-study cardiac record created in this work. For each ECG study, we collect the

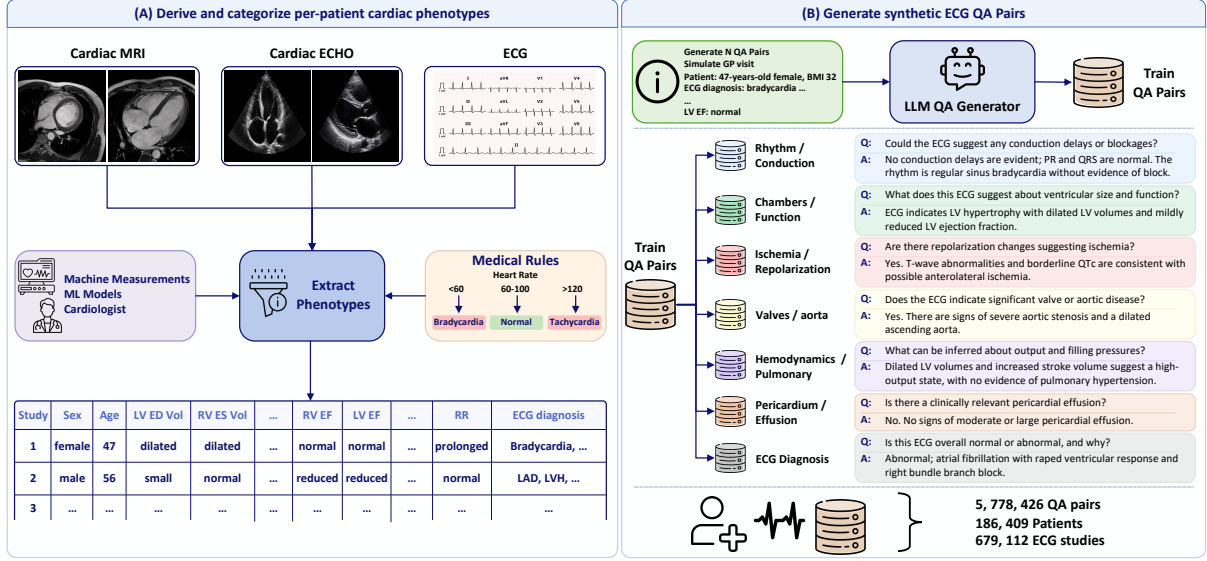


Figure 10: Clinical QA supervision from cardiac records. (A) ECG, CMR, echocardiography, demographic, and symptom information are converted into structured per-study cardiac records. (B) A locally hosted instruction-tuned LLM uses each structured record to generate ECG-focused QA pairs that simulate front-line clinical questions.

available information linked to that ECG, including ECG measurements, diagnostic statements or reports, imaging-derived phenotypes when available, and patient demographics or symptoms when available. These data are converted into a structured text record. A QA generator then uses the record to produce clinically grounded front-line cardiac questions and answers. The generator does not see the raw ECG waveform. After generation, each QA pair is paired with the corresponding ECG waveform for ECG-LLM training, so that ECG-LLM learns to answer from ECG tokens and question tokens at inference time.

The structured records combine complementary information from the four cohorts. MIMIC-IV-ECG [11] provides device-derived ECG measurements and diagnostic statements, while PTB-XL [12] provides cardiologist ECG reports and structured SCP-ECG diagnostic labels. UK Biobank [10] provides device-derived ECG measurements together with CMR-derived cardiac phenotypes extracted from long-axis and short-axis cardiac MRI using a machine-learning segmentation pipeline [65]. EchoNext [13] provides device-derived ECG measurements together with echocardiographic measurements and binary structural heart disease labels. The complete list of dataset-specific ECG, CMR, echocardiographic, demographic, symptom, and report-derived fields used to construct the structured records is provided in Supplementary Table 13.

Before prompting the QA generator, continuous measurements are converted into clinically meaningful categorical descriptors. This makes the records closer to clinical language, reduces dependence on exact numerical values, and avoids asking the generator to perform thresholding inside the prompt. Prior work has shown that language models can be unreliable on exact numerical reasoning compared with symbolic or textual representations [66]. We therefore apply rule-based categorization across all four cohorts using clinical ECG thresholds, sex-specific CMR reference ranges, ventricular-function thresholds, wall-thickness categories, strain reference ranges, and age- and sex-specific aortic-distensibility ranges. When sex, age, or body surface area is available, the corresponding sex-specific, age-specific, or indexed reference rule is applied. The complete rule set is provided in Supplementary Table 14.

Each structured record is prompted independently using GPT-OSS-20B [50, 51, 52]. The prompt describes a front-line care scenario in which a clinician asks clinically meaningful ECG-based questions and an expert ECG interpreter answers using the cardiac record. The generator

is instructed to ask ECG-focused questions grounded in the record, express imaging-derived descriptors as ECG-inferable cardiac findings without explicitly mentioning CMR, echocardiography, or imaging, and produce answers faithful to the supplied information. This design lets imaging-derived supervision teach the ECG-conditioned model about broader cardiac phenotypes while keeping inference ECG-only. To encourage lexical diversity, we randomize the order of phenotype fields and add a study-specific randomization token before prompting the QA generator. The generation procedure is applied independently to each ECG study and yields triplets (x, q, a) , where x is the ECG waveform, q is a front-line cardiac question, and a is the corresponding grounded free-text answer.

For UK Biobank and EchoNext, where ECGs are paired with imaging-derived phenotypes, we also use the same structured-record setup to generate chain-of-thought QA pairs. These examples are intended to encourage clinically plausible intermediate reasoning between ECG-visible findings and downstream structural or functional cardiac phenotypes, while keeping the final answer grounded in the same record. Prior work has shown that intermediate reasoning traces can improve complex reasoning in large language models and can serve as useful supervision during rationale-based distillation [67, 68]. For the UK Biobank specifically, where most participants are normal for many phenotypes, we additionally generate an abnormal-enriched subset. We retain only studies with at least one abnormal imaging-derived cardiac phenotype. For each retained study, we create an abnormal-focused record by keeping only non-normal findings, including abnormal CMR phenotype descriptors, abnormal ECG measurements, and ECG diagnosis text when available. This increases exposure to clinically important abnormal findings and partially mitigates the dominance of normal cases in UK Biobank.

In addition to structured-record QA, we include two external ECG-language sources to strengthen conventional ECG interpretation. First, we use ECG-QA training examples [37], which provide established ECG-conditioned question-answering tasks with candidate options. These examples anchor the model to standard ECG measurement and diagnostic QA. Second, we use selected items from the ECGInstruct dataset [47]. We retain open-ended and close-ended questions from feature, rhythm, morphology, and report tasks because they correspond to clinically interpretable ECG findings, rhythm assessment, waveform morphology, and summary reporting. These tasks complement the structured-record QA without shifting the model away from front-line cardiac question answering. ECG-Instruct prompts are lightly cleaned to remove image placeholders and role-play preambles and to replace image-specific wording with ECG-specific wording. All external QA examples are converted to the same ECG-LLM chat format used for structured-record QA (Supplementary Section 11). The final instruction-tuning corpus composition by dataset, QA source, and QA subset or task is provided in Supplementary Table 6.

4.5 Training and implementation details

We train ECG-LLM on the ECG-specific QA corpus described above. The ECG encoder described in Section 4.2 is first pretrained as a ViT-MAE [59] on ECGs from the four training cohorts: UK Biobank [10], MIMIC-IV-ECG [11], PTB-XL [12], and EchoNext [13]. During MAE pretraining, ECGs are randomly cropped to 2,500 time steps and augmented with Fourier-transform surrogate perturbation, jitter, and rescaling. The augmented ECGs are tokenized using patches of size 1×100 , yielding 300 ECG tokens before masking. The MAE objective masks 75% of tokens and reconstructs the masked patches. The pretraining batch size is 768, and optimisation uses bfloat16 mixed precision. Full ECG encoder pretraining details are provided in Supplementary Table 15.

After ECG encoder pretraining, we train ECG-LLM for ECG-conditioned question answering. The pretrained ECG encoder and original LLM backbone remain frozen. The trainable parameters are the ECG projector and LoRA adapter weights [63]. LoRA adapters are attached to the query, key, value, and output projections of self-attention layers, with rank 128, scaling factor 128, and dropout 0.1. We use a two-phase curriculum to stabilise ECG-language

alignment. The first phase lets the projector move ECG tokens into an approximate region of the LLM embedding space before LLM-side adaptation begins. This reduces the chance that early joint updates force the LoRA adapters to compensate for poorly scaled ECG tokens or rely mainly on the LLM’s text priors instead of learning ECG-conditioned representations. For the first 1,500 optimiser steps, only the projector is trained, after which the projector and LoRA adapters are trained jointly under the causal language-modelling objective. We optimise trainable parameters with paged AdamW 8-bit [69]. The projector learning rate is 10^{-4} , the LoRA learning rate is 10^{-5} , and weight decay is 10^{-3} . We use linear warm-up over 3% of optimiser steps followed by cosine decay to 0.3 of the peak learning rate. ECG-LLM training is performed for one epoch using bfloat16 precision on four NVIDIA A100-SXM4-80GB GPUs. The per-device batch size is 8 with gradient accumulation over 4 steps, giving an effective batch size of 128 ECG–question–answer examples. During ECG-LLM training, ECGs are randomly cropped to 2,500 time steps, and no additional ECG augmentations are applied. The maximum tokenized sequence length is 1,200 tokens, covering 99.5% of training examples. Longer examples are truncated. Full implementation details are provided in Supplementary Table 16.

Data availability

The source datasets analysed in this study are available through public or controlled-access repositories. UK Biobank data are available to approved researchers following application through the UK Biobank Research Analysis Platform (link). For this study, permission to access and analyze the UK Biobank data was approved under application 87802, with initial approval granted in July 2022. MIMIC-IV-ECG is available to credentialed users through PhysioNet under the applicable data-use agreement (link). PTB-XL is publicly available through PhysioNet (link), and EchoNext is available through PhysioNet under its applicable access conditions (link). The original source data cannot be redistributed by the authors because they remain subject to the access conditions and data-use agreements of the respective data providers. Derived ECG question–answer supervision data, cohort-specific metadata mappings, and evaluation splits generated in this study are planned for release upon publication to the extent permitted by the corresponding data-use agreements.

Code availability

The underlying code for this study, including preprocessing, question–answer generation, model training and fine-tuning, inference, evaluation, and reproducibility scripts, will be available in the GitHub repository at <https://github.com/placeholder/> upon publication. Model checkpoints for the main experiments will also be made available through the repository, subject to the applicable licences and data-use agreements.

Acknowledgments

This research has been conducted using the UK Biobank Resource under Application Number 87802. This work is funded by the European Research Council (ERC) project Deep4MI (884622).

Author contributions

A.S. and D.R. conceived and designed the study. A.S. developed the methodology, implemented the models, performed the experiments and analyses, and wrote the original manuscript. F.J., J.K., K.-L.L., and E.M. provided clinical expertise and contributed to the interpretation of the

results. D.R. supervised the study and provided technical and scientific guidance. All authors critically reviewed, revised, and approved the final manuscript.

Competing interests

The authors declare no competing interests.

References

- [1] World Health Organization. WHO Cardiovascular-Diseases, 2025.
- [2] Global Burden of Cardiovascular Diseases and Risks 2023 Collaborators. Global, regional, and national burden of cardiovascular diseases and risk factors in 204 countries and territories, 1990-2023. *Journal of the American College of Cardiology*, 86(22):2167–2243, 2025.
- [3] R. Sacha Bhatia et al. Electrocardiograms in low-risk patients undergoing an annual health examination. *JAMA Internal Medicine*, 177(9):1326–1333, Sep 2017.
- [4] Quirino Ciampi and Bruno Villari. Role of echocardiography in diagnosis and risk stratification in heart failure with left ventricular systolic dysfunction. *Cardiovascular Ultrasound*, 5(1):34, Oct 2007.
- [5] G Habib and A Torbicki. The role of echocardiography in the diagnosis and management of patients with pulmonary hypertension. *Eur. Respir. Rev.*, 19(118):288–299, Dec 2010.
- [6] North J Noelck, Heather A Perry, Phyllis L Talley, and D Elizabeth Le. A narrative review of the clinical applications of echocardiography in right heart failure. *J. Clin. Med.*, 14(15):5505, Aug 2025.
- [7] Michael Salerno et al. Recent advances in cardiovascular magnetic resonance: Techniques and applications. *Circ. Cardiovasc. Imaging*, 10(6), Jun 2017.
- [8] Henry Anselmo Mayala, Khamis Hassan Bakari, and Zhaohui Wang. The role of Cardiac Magnetic Resonance (CMR) in the diagnosis of cardiomyopathy: A systematic review. *Malawi Med. J.*, 31(3):241–245, Sep 2019.
- [9] Waseem Hindieh, Raymond Chan, and Harry Rakowski. Complementary role of echocardiography and cardiac magnetic resonance in hypertrophic cardiomyopathy. *Curr. Cardiol. Rep.*, 19(9):81, Sep 2017.
- [10] C. Sudlow et al. UK Biobank: An open access resource for identifying the causes of a wide range of complex diseases of middle and old age. *PLoS Med.*, 12(3):e1001779, 2015.
- [11] Brian Gow et al. MIMIC-IV-ECG: Diagnostic Electrocardiogram Matched Subset. *PhysioNet*, Sep 2023. Version 1.0.
- [12] Patrick Wagner, Nils Strodthoff, Ralf-Dieter Bousseljot, Wojciech Samek, and Tobias Schaeffer. PTB-XL, a large publicly available electrocardiography dataset. *PhysioNet*, Nov 2022. Version 1.0.3.
- [13] Pierre Elias and Joshua Finer. EchoNext: A Dataset for Detecting Echocardiogram-Confirmed Structural Heart Disease from ECGs. *PhysioNet*, Sep 2025. Version 1.1.0.
- [14] Awni Y. Hannun et al. Cardiologist-level arrhythmia detection and classification in ambulatory electrocardiograms using a deep neural network. *Nature Medicine*, 25(1):65–69, Jan 2019.

- [15] Antônio H. Ribeiro et al. Automatic diagnosis of the 12-lead ECG using a deep neural network. *Nature Communications*, 11(1):1760, Apr 2020.
- [16] Chenyang Qiu, Hao Li, Chaoqun Qi, and Bo Li. Enhancing ECG classification with continuous wavelet transform and multi-branch transformer. *Heliyon*, 10(5):e26147, 2024.
- [17] Zachi I. Attia et al. Screening for cardiac contractile dysfunction using an artificial intelligence-enabled electrocardiogram. *Nature Medicine*, 25(1):70–74, Jan 2019.
- [18] Sushravya Raghunath et al. Deep Neural Networks Can Predict New-Onset Atrial Fibrillation From the 12-Lead ECG and Help Identify Those at Risk of Atrial Fibrillation-Related Stroke. *Circulation*, 143(13):1287–1298, 2021.
- [19] Michal Cohen-Shelly et al. Electrocardiogram screening for aortic valve stenosis using artificial intelligence. *Eur. Heart J.*, 42(30):2885–2896, Aug 2021.
- [20] Demilade Adedinsewo et al. Artificial intelligence-enabled ECG algorithm to identify patients with left ventricular systolic dysfunction presenting to the emergency department with dyspnea. *Circ. Arrhythm. Electrophysiol.*, 13(8):e008437, Aug 2020.
- [21] Wei-Yin Ko et al. Detection of hypertrophic cardiomyopathy using a convolutional neural network-enabled electrocardiogram. *J. Am. Coll. Cardiol.*, 75(7):722–733, Feb 2020.
- [22] Geoffrey H. Tison et al. Assessment of disease status and treatment response with artificial intelligence-enhanced electrocardiography in obstructive hypertrophic cardiomyopathy. *J. Am. Coll. Cardiol.*, 79(10):1032–1034, Mar 2022.
- [23] Eunjung Lee et al. Artificial intelligence-enabled ECG for left ventricular diastolic function and filling pressure. *npj Digital Medicine*, 7(1):4, Jan 2024.
- [24] Sam F. Friedman et al. Unsupervised deep learning of electrocardiograms enables scalable human disease profiling. *npj Digital Medicine*, 8(1):23, Jan 2025.
- [25] Mostafa A. Al-Alusi et al. A deep learning digital biomarker to detect hypertension and stratify cardiovascular risk from the electrocardiogram. *npj Digital Medicine*, 8(1):120, Feb 2025.
- [26] Alexander Selivanov, Philip Müller, Özgün Turgut, Nil Stolt-Ansó, and Daniel Rueckert. Global and Local Contrastive Learning for Joint Representations from Cardiac MRI and ECG . In *proceedings of Medical Image Computing and Computer Assisted Intervention – MICCAI 2025*, volume 15960 of *Lecture Notes in Computer Science*, pages 217–227. Springer Nature Switzerland, Sep 2025.
- [27] Yuan Gao, Sangwook Kim, and Chris McIntosh. EchoingECG: An Electrocardiogram Cross-Modal Model for Echocardiogram Tasks . In *proceedings of Medical Image Computing and Computer Assisted Intervention – MICCAI 2025*, volume 15960 of *Lecture Notes in Computer Science*, pages 175–185. Springer Nature Switzerland, Sep 2025.
- [28] Robbie Holland et al. Specialized curricula for training vision language models in retinal image analysis. *npj Digital Medicine*, 8(1):532, Aug 2025.
- [29] Zhihong Chen et al. Chexagent: Towards a foundation model for chest x-ray interpretation. *arXiv preprint arXiv:2401.12208*, 2024.
- [30] Chunyuan Li et al. LLaVA-Med: Training a Large Language-and-Vision Assistant for Biomedicine in One Day. In *Advances in Neural Information Processing Systems*, volume 36, pages 28541–28564, 2023.

- [31] Michael Moor et al. Med-Flamingo: a Multimodal Medical Few-shot Learner. In Stefan Hegselmann, Antonio Parziale, Divya Shanmugam, Shengpu Tang, Mercy Nyamewaa Asiedu, Serina Chang, Tom Hartvigsen, and Harvineet Singh, editors, *Proceedings of the 3rd Machine Learning for Health Symposium*, volume 225 of *Proceedings of Machine Learning Research*, pages 353–367. PMLR, Dec 2023.
- [32] Xiao Gu et al. Cardiac health assessment across scenarios and devices using a multimodal foundation model pretrained on data from 1.7 million individuals. *Nature Machine Intelligence*, 8(2):220–233, Feb 2026.
- [33] Yubao Zhao, Jiaju Kang, Tian Zhang, Puyu Han, and Tong Chen. ECG-Chat: A Large ECG-Language Model for Cardiac Disease Diagnosis. In *2025 IEEE International Conference on Multimedia and Expo (ICME)*, pages 1–6, 2025.
- [34] Ruoqi Liu, Yuelin Bai, Xiang Yue, and Ping Zhang. Teaching multimodal LLMs to comprehend 12-lead electrocardiographic images. *npj Digital Medicine*, 9(1):349, Mar 2026.
- [35] Dongsheng Tian et al. ECG-Doctor: An Interpretable Multimodal ECG Diagnosis Framework Based on Large Language Models. In *Proceedings of the 34th ACM International Conference on Information and Knowledge Management, CIKM '25*, pages 2863–2873, New York, NY, USA, 2025. Association for Computing Machinery.
- [36] Xiang Lan et al. GEM: Empowering MLLM for Grounded ECG Understanding with Time Series and Images. In *The Thirty-ninth Annual Conference on Neural Information Processing Systems*, 2025.
- [37] Jungwoo Oh, Gyubok Lee, Seongsu Bae, Joon-myung Kwon, and Edward Choi. ECG-QA: A comprehensive question answering dataset combined with electrocardiogram. *Advances in Neural Information Processing Systems*, 36:66277–66288, 2023.
- [38] Jong Hak Moon, Hyungyung Lee, Woncheol Shin, Young-Hak Kim, and Edward Choi. Multi-modal understanding and generation for medical images and text via vision-language pre-training. *IEEE Journal of Biomedical and Health Informatics*, 26(12):6070–6080, 2022.
- [39] Zhihong Chen et al. Multi-modal Masked Autoencoders for Medical Vision-and-Language Pre-training. In *Medical Image Computing and Computer Assisted Intervention – MICCAI 2022*, volume 13435 of *Lecture Notes in Computer Science*, pages 679–689, Cham, 2022. Springer Nature Switzerland.
- [40] Andrew B. Sellergren et al. MedGemma Technical Report, 2025.
- [41] Google DeepMind. MedGemma. <https://deepmind.google/models/gemma/medgemma/>, 2025. Model overview.
- [42] Google. google/medgemma-4b-it. <https://huggingface.co/google/medgemma-4b-it>, 2025. Hugging Face model card.
- [43] Google Health AI Developer Foundations. MedGemma. <https://developers.google.com/health-ai-developer-foundations/medgemma>, 2025. Developer documentation.
- [44] PULSE-ECG. PULSE-7B. <https://huggingface.co/PULSE-ECG/PULSE-7B>, 2024. Hugging Face model card.
- [45] Darren Edge et al. From local to global: A graph RAG approach to query-focused summarization, 2024.

- [46] Dilara Soyly, Christopher Potts, and Omar Khattab. Fine-tuning and prompt optimization: Two great steps that work better together. In Yaser Al-Onaizan, Mohit Bansal, and Yun-Nung Chen, editors, *Proceedings of the 2024 Conference on Empirical Methods in Natural Language Processing*, pages 10696–10710, Miami, Florida, USA, Nov 2024. Association for Computational Linguistics.
- [47] PULSE-ECG. ECGInstruct. <https://huggingface.co/datasets/PULSE-ECG/ECGInstruct>, 2024. Hugging Face dataset card.
- [48] Kshama Kodthalu Shivashankara et al. ECG-Image-Kit: a synthetic image generation toolbox to facilitate deep learning-based electrocardiogram digitization. *Physiological Measurement*, 45(5):055019, May 2024.
- [49] AlphaNumerics Lab. ECG-Image-Kit: A toolkit for synthesis, analysis, and digitization of electrocardiogram images. <https://github.com/alphanumericslab/ecg-image-kit>, 2024. GitHub repository.
- [50] OpenAI. Introducing gpt-oss. <https://openai.com/index/introducing-gpt-oss/>, 2025.
- [51] OpenAI. gpt-oss-20b model card. <https://huggingface.co/openai/gpt-oss-20b>, 2025.
- [52] OpenAI. gpt-oss-120b and gpt-oss-20b model card, 2025.
- [53] Lianmin Zheng et al. Judging LLM-as-a-Judge with MT-Bench and Chatbot Arena. In *Advances in Neural Information Processing Systems*, volume 36, pages 46595–46623, 2023.
- [54] Zilong Wang, Xufang Luo, Xinyang Jiang, Dongsheng Li, and Lili Qiu. LLM-RadJudge: Achieving Radiologist-Level Evaluation for X-Ray Report Generation, 2024.
- [55] Paul Hager et al. Evaluation and mitigation of the limitations of large language models in clinical decision-making. *Nature Medicine*, 30(9):2613–2622, Sep 2024.
- [56] Deyao Zhu, Jun Chen, Xiaoqian Shen, Xiang Li, and Mohamed Elhoseiny. MiniGPT-4: Enhancing Vision-Language Understanding with Advanced Large Language Models. In *The Twelfth International Conference on Learning Representations*, 2024.
- [57] Haotian Liu, Chunyuan Li, Qingyang Wu, and Yong Jae Lee. Visual Instruction Tuning. In A. Oh, T. Naumann, A. Globerson, K. Saenko, M. Hardt, and S. Levine, editors, *Advances in Neural Information Processing Systems*, volume 36, pages 34892–34916, 2023.
- [58] Haotian Liu, Chunyuan Li, Yuheng Li, and Yong Jae Lee. Improved Baselines with Visual Instruction Tuning. In *2024 IEEE/CVF Conference on Computer Vision and Pattern Recognition (CVPR)*, pages 26286–26296, 2024.
- [59] Kaiming He et al. Masked Autoencoders Are Scalable Vision Learners. In *2022 IEEE/CVF Conference on Computer Vision and Pattern Recognition (CVPR)*, pages 15979–15988. IEEE, Jun 2022.
- [60] Ashish Vaswani et al. Attention is All you Need. In I. Guyon, U. Von Luxburg, S. Bengio, H. Wallach, R. Fergus, S. Vishwanathan, and R. Garnett, editors, *Advances in Neural Information Processing Systems*, volume 30, pages 5998–6008. Curran Associates, Inc., 2017.
- [61] AI at Meta. Introducing Llama 3.1: Our most capable models to date. <https://ai.meta.com/blog/meta-llama-3-1/>, 2024.

- [62] Meta AI. Llama 3.1 8B Instruct model card. <https://huggingface.co/meta-llama/Llama-3.1-8B-Instruct>, 2024. Model card.
- [63] Edward J. Hu et al. LoRA: Low-Rank Adaptation of Large Language Models. In *International Conference on Learning Representations*, 2022.
- [64] David Cuesta Frau, Daniel Novák, Vladimir Eck, Juan C Pérez Cortés, and Gabriela Andreu García. Electrocardiogram baseline removal using wavelet approximations. In *Proceeding of the 15th Biennial Eurasip Conference BIOSIGNAL*, pages 136–138, 2000.
- [65] Wenjia Bai et al. A population-based phenome-wide association study of cardiac and aortic structure and function. *Nat. Med.*, 26(10):1654–1662, Oct 2020.
- [66] Gilad Yehudai, Haim Kaplan, Asma Ghandeharioun, Mor Geva, and Amir Globerson. When Can Transformers Count to n ?, 2024.
- [67] Jason Wei et al. Chain-of-Thought Prompting Elicits Reasoning in Large Language Models. In *Advances in Neural Information Processing Systems*, volume 35, pages 24824–24837, 2022.
- [68] Cheng-Yu Hsieh et al. Distilling Step-by-Step! Outperforming Larger Language Models with Less Training Data and Smaller Model Sizes. In Anna Rogers, Jordan Boyd-Graber, and Naoaki Okazaki, editors, *Findings of the Association for Computational Linguistics: ACL 2023*, pages 8003–8017, Toronto, Canada, Jul 2023. Association for Computational Linguistics.
- [69] Tim Dettmers, Mike Lewis, Sam Shleifer, and Luke Zettlemoyer. 8-bit Optimizers via Block-wise Quantization. In *International Conference on Learning Representations*, 2022.
- [70] Akshay Goel. Langextract v1.4.0, 2026. Version 1.4.0.
- [71] Google. Langextract. <https://github.com/google/langextract>, 2025. GitHub repository.

Supplementary Information

1 Dataset statistics

Table 5 shows the number of ECG recordings and patients in each dataset split. Table 6 summarizes the final ECG-LLM instruction-tuning corpus.

Split	Entity	UK Biobank	MIMIC-IV-ECG	PTB-XL	EchoNext	All
Train	ECG studies	34,280	560,791	11,566	72,475	679,112
Train	Patients	34,280	115,508	10,403	26,218	186,409
Validation	ECG studies	–	62,251	1,274	4,626	68,151
Validation	Patients	–	12,834	1,155	4,626	18,615
Test	ECG studies	6,049	156,851	3,214	5,442	171,556
Test	Patients	6,049	32,086	2,891	5,442	46,468

Table 5: Dataset split statistics. The table reports the number of ECG studies and unique patients in the training, validation, and test splits for each ECG cohort. UK Biobank was split subject-wise into training and test sets only.

Dataset	QA source	Type	# Samples
UK Biobank	Structured-record QA	ECG and CMR, standard	670,545
	Structured-record QA	ECG and CMR, abnormal-enriched	173,238
	Structured-record QA	ECG and CMR, chain-of-thought	167,701
EchoNext	Structured-record QA	ECG and echo, standard	546,156
	Structured-record QA	ECG and echo, chain-of-thought	129,625
MIMIC-IV-ECG	Structured-record QA	ECG interpretation, standard	2,775,751
	ECG-QA	Single-Verify	125,380
	ECG-QA	Single-Choose	107,860
	ECG-QA	Single-Query	153,156
	ECGInstruct	Feature (Close/Open)	3,865
	ECGInstruct	Rhythm (Close/Open)	38,124
	ECGInstruct	Morphology (Close/Open)	62,663
	ECGInstruct	Report (Open)	345,561
PTB-XL	Structured-record QA	ECG interpretation, standard	127,219
	ECG-QA	Single-Verify	62,554
	ECG-QA	Single-Choose	50,015
	ECG-QA	Single-Query	46,737
	ECGInstruct	Feature (Close/Open)	57,108
	ECGInstruct	Rhythm (Close/Open)	18,899
	ECGInstruct	Morphology (Close/Open)	108,379
	ECGInstruct	Report (Open)	7,890
Total			5,778,426

Table 6: Final ECG-LLM instruction-tuning corpus. Structured-record QA denotes generated examples produced from the per-study cardiac records created in this work.

2 UK Biobank ECG and CMR phenotype question bank

Question group	ECG-LLM questions	Target phenotypes
Heart rate and rhythm	What does the ECG reveal about the patient's heart rate and rhythm? Review the ECG and describe the rhythm and rate.	ECG heart rate [bpm]; RR interval [ms].
P-wave characteristics and atrioventricular intervals	Review this ECG and describe atrial activity and AV conduction. Review the ECG and infer the P-wave characteristics and atrioventricular intervals.	P duration [ms]; PP interval [ms]; PQ interval [ms].
QRS duration and ventricular conduction	Is the QRS duration within normal limits? Review the ECG and identify any conduction delays.	QRS duration [ms].
QT interval and ventricular repolarization	Are there any repolarization abnormalities on this ECG? Is the QT interval within normal limits?	QT interval [ms]; QTC interval [ms].
Electrical axes	Are the electrical axes within normal limits? Does the ECG show any axis deviation?	P axis [degrees]; R axis [degrees]; T axis [degrees].
Left ventricular systolic function	Review the ECG and describe LV ejection fraction. What does this ECG reveal about left ventricular systolic function?	LV ejection fraction.
Left ventricular stroke volume	Review the ECG and describe LV stroke volume. What does this ECG reveal about left ventricular stroke volume?	LV stroke volume.
Left ventricular volumes	Review the ECG and describe any abnormalities found for LV. What does this ECG reveal about left ventricular volumes?	LV end diastolic volume; LV end systolic volume.
Left ventricular myocardial mass and wall thickness	Is there any evidence of ventricular hypertrophy or increased myocardial mass? What does this ECG reveal about left ventricular wall thickness?	LV myocardial mass; LV mean myocardial wall thickness global.
Left ventricular strain	Do the strain measurements suggest any myocardial dysfunction? Review this ECG and describe any abnormalities in left ventricular strain.	LV circumferential strain global; LV longitudinal strain global; LV radial strain global.
Cardiac output and cardiac index	What does this ECG indicate about cardiac output and cardiac index? Examine this ECG and tell if the cardiac output or index within normal limits?	LV cardiac output; cardiac index [litres/min/m2].
Right ventricular stroke volume	Review this ECG and describe RV stroke volume. What does this ECG reveal about right ventricular stroke volume?	RV stroke volume.
Right ventricular ejection fraction	Review this ECG and describe RV ejection fraction. What does this ECG reveal about right ventricular systolic function?	RV ejection fraction.
Right ventricular volumes and function	Review this ECG and describe all abnormalities for RV. Review the ECG and indicate about abnormalities in right ventricular volumes.	RV end diastolic volume; RV end systolic volume.
Left atrial volumes and function	What does this ECG suggest about left atrial volumes and function? Examine this ECG and tell if there is any evidence of left atrial enlargement or dysfunction?	LA maximum volume; LA minimum volume; LA stroke volume; LA ejection fraction.
Right atrial volumes and function	What does this ECG reveal about right atrial volumes and function? Examine this ECG and tell if there is any evidence of right atrial enlargement or dysfunction?	RA maximum volume; RA minimum volume; RA stroke volume; RA ejection fraction.
Aortic distensibility	Review this ECG and state whether there is any indirect evidence that could suggest abnormal aortic or vascular stiffness.	Ascending aorta distensibility; Descending aorta distensibility.

Table 7: UK Biobank question bank used for ECG-derived and CMR-derived phenotype elicitation from ECG alone.

3 UK Biobank phenotype extraction protocol

Because ECG-LLM was trained to generate free-text answers, the UK Biobank phenotype evaluation required an additional mapping step from generated text to the predefined phenotype-category schema. We used the Google LangExtract framework [70, 71] to apply a schema-constrained extraction prompt to ECG-LLM answers. LangExtract is a Python library for extracting structured information from unstructured text using LLMs and user-defined extraction instructions. In our setup, the backend model was a locally deployed GPT-OSS-20B instance [50, 51, 52].

For each ECG and question group, the extractor received the generated answer, the target phenotype fields assigned to that question group, and the allowed category labels for each target phenotype. It was instructed to return allowed category labels only for those target fields found in the ECG-LLM answers. This group-specific extraction kept the evaluation aligned with the clinical intent of the question: the model was neither rewarded nor penalized for mentioning unrelated phenotypes outside the asked group. This step was necessary because a free-text answer could describe several related phenotypes together, refer to multiple measurements collectively, or express a category without repeating the exact phenotype name used in the evaluation schema. Therefore, simple string matching or regular expressions would not reliably map answers to field-specific category labels. The extractor served only as a mapping layer from free text to the predefined category space used for evaluation.

Prompt for schema-constrained phenotype extraction

```
Extract discrete cardiac phenotype categories.
Use only the phenotype field names listed below as possible extraction classes.
Use the allowed category words as the extraction text.
Do not extract numeric values or explanations.
If the answer says several named phenotypes are normal, extract normal for each named
↪ phenotype.
Sometimes, phenotype fields may not contain dimension information in square brackets, this is
↪ still valid.
If the answer says one phenotype is abnormal but does not give a discrete allowed category,
↪ do not extract it.
ALLOWED PHENOTYPES AND CATEGORIES:
- <phenotype field 1>: <allowed category 1>, <allowed category 2>, ...
- <phenotype field 2>: <allowed category 1>, <allowed category 2>, ...
- ...
```

4 UK Biobank phenotype extraction results

Field	n	K	Cov.	Miss.	Acc.	BAcc.	wF1	Maj.	mF1
<i>ECG measurements</i>									
ECG heart rate [bpm]	1000	5	1.000	0.000	0.958	0.974	0.958	0.126	0.977
R axis [degrees]	1000	3	1.000	0.000	0.987	0.919	0.987	0.315	0.940
RR interval [ms]	1000	5	1.000	0.000	0.944	0.969	0.944	0.127	0.926
P duration [ms]	979	2	1.000	0.000	0.974	0.824	0.972	0.481	0.886
QRS duration [ms]	1000	3	0.987	0.013	0.976	0.886	0.976	0.320	0.874
T axis [degrees]	1000	4	1.000	0.000	0.888	0.681	0.880	0.212	0.726
QTC interval [ms]	1000	3	1.000	0.000	0.962	0.652	0.958	0.324	0.719
PP interval [ms]	999	5	1.000	0.000	0.930	0.605	0.927	0.126	0.637
P axis [degrees]	938	3	1.000	0.000	0.908	0.540	0.888	0.310	0.614
PQ interval [ms]	980	3	1.000	0.000	0.906	0.458	0.879	0.313	0.522
<i>CMR-derived phenotypes</i>									
LV mean myocardial wall thickness global	990	2	1.000	0.000	0.965	0.974	0.966	0.445	0.947
RV end diastolic volume	1000	3	0.999	0.001	0.993	0.708	0.991	0.327	0.708
LV radial strain global	833	3	0.975	0.025	0.998	0.667	0.996	0.311	0.666
RV end systolic volume	1000	3	0.999	0.001	0.995	0.667	0.993	0.327	0.641
RA minimum volume	988	3	1.000	0.000	0.910	0.594	0.885	0.314	0.592
RA maximum volume	989	3	0.999	0.001	0.906	0.569	0.894	0.312	0.560
LA maximum volume	989	3	0.999	0.001	0.950	0.473	0.937	0.322	0.547
LV end systolic volume	998	3	1.000	0.000	0.943	0.460	0.930	0.318	0.507
LV myocardial mass	1000	2	0.977	0.023	0.950	0.495	0.934	0.488	0.487
Ascending aorta distensibility	872	2	0.814	0.186	0.866	0.484	0.891	0.489	0.486
RV ejection fraction	999	5	1.000	0.000	0.981	0.443	0.974	0.197	0.469
LA minimum volume	988	3	1.000	0.000	0.881	0.429	0.850	0.305	0.459
RV stroke volume	999	3	1.000	0.000	0.961	0.472	0.955	0.327	0.458
RA stroke volume	989	3	0.999	0.001	0.884	0.403	0.866	0.314	0.416
LV end diastolic volume	1000	3	0.998	0.002	0.945	0.383	0.925	0.322	0.410
Cardiac index [litres/min/m ²]	843	5	1.000	0.000	0.696	0.396	0.672	0.114	0.397
LA ejection fraction	989	3	0.999	0.001	0.908	0.369	0.868	0.316	0.381
LV ejection fraction	1000	5	1.000	0.000	0.946	0.321	0.934	0.191	0.371
LV stroke volume	1000	3	1.000	0.000	0.941	0.360	0.917	0.323	0.371
RA ejection fraction	989	3	0.999	0.001	0.876	0.361	0.831	0.310	0.363
LV longitudinal strain global	815	3	0.996	0.004	0.966	0.346	0.950	0.327	0.352
Descending aorta distensibility	872	3	0.810	0.190	0.857	0.358	0.887	0.326	0.339
LA stroke volume	989	3	0.999	0.001	0.956	0.333	0.935	0.326	0.326
LV circumferential strain global	833	3	0.975	0.025	0.950	0.333	0.925	0.324	0.325
LV cardiac output	1000	5	0.843	0.157	0.766	0.271	0.695	0.170	0.277

Table 8: Full UK Biobank phenotype extraction results for ECG-derived and CMR-derived fields. The table reports field-level sample size (n), number of possible classification classes (K), coverage, missing rate, accuracy, balanced accuracy, weighted-F1, majority-class macro-F1, and macro-F1. Fields are ordered by macro-F1 within each group.

5 EchoNext ECG-only binary phenotype questions

Phenotype	ECG-LLM question
$LVEF \leq 45\%$	Does this ECG suggest moderately or severely reduced left ventricular systolic function?
$LVWT \geq 1.3$ cm	Does this ECG show moderate left ventricular hypertrophy with the maximum of the interventricular septum (IVS) or posterior wall (LVPW) thickness greater than or equal to 1.3 cm?
$AS \geq$ moderate	Does this ECG show moderate or severe aortic stenosis?
$AR \geq$ moderate	Does this ECG suggest moderate or severe aortic regurgitation?
$MR \geq$ moderate	Does this ECG suggest moderate or severe mitral regurgitation?
$TR \geq$ moderate	Does this ECG suggest moderate or severe tricuspid regurgitation?
$PR \geq$ moderate	Does this ECG suggest moderate or severe pulmonary regurgitation?
RV systolic dysf. $\geq$ moderate	Does this ECG suggest moderate or severe right ventricular systolic dysfunction?
Pericardial eff. $\geq$ moderate/large	Does this ECG suggest a moderate or large pericardial effusion?
$PASP \geq 45$ mmHg	Does this ECG suggest pulmonary hypertension, as would be expected with a pulmonary artery systolic pressure of at least 45 mmHg?
TR Vmax ≥ 3.2 m/s	Does this ECG suggest a tricuspid regurgitation jet velocity of at least 3.2 m/s?

Table 9: EchoNext phenotype-specific questions used for ECG-only yes/no question answering. For each ECG and phenotype, ECG-LLM first answered the initial question in free text. In the same conversation, the model was then asked: “Give the final answer to the original question as yes or no.” The final yes/no response was parsed and compared with the corresponding echocardiographic binary label.

6 PTB-XL ECG-QA performance by question format and attribute

Question format	Attribute type	n	EM	μP	μR	$\mu F1$
Single-Verify	SCP code	8,243	79.52	79.52	79.52	79.52
Single-Verify	Heart axis	237	89.87	89.87	89.87	89.87
Single-Verify	Stage of infarction	168	82.74	82.74	82.74	82.74
Single-Verify	Noise	3,183	67.73	67.73	67.73	67.73
Single-Verify	Extra systole	230	77.39	77.39	77.39	77.39
Single-Verify	Numeric feature	1,020	77.84	77.84	77.84	77.84
Single-Choose	SCP code	9,358	69.30	77.50	72.95	75.16
Single-Choose	Heart axis	36	75.00	75.00	75.00	75.00
Single-Choose	Stage of infarction	18	55.56	55.56	55.56	55.56
Single-Choose	Noise	317	43.22	55.21	46.27	50.35
Single-Choose	Extra systole	22	50.00	54.55	46.15	50.00
Single-Choose	Numeric feature	104	60.58	60.58	60.58	60.58
Single-Query	SCP code	9,208	43.07	59.47	57.99	58.72
Single-Query	Heart axis	167	74.85	74.85	74.85	74.85
Single-Query	Stage of infarction	178	46.07	46.07	46.07	46.07
Single-Query	Noise	3,960	38.86	44.48	46.04	45.24
Single-Query	Extra systole	241	54.77	56.85	55.02	55.92
Single-Query	Numeric feature	4,403	25.03	65.75	65.96	65.85

Table 10: PTB-XL ECG-QA performance stratified by question format and ECG-QA attribute type. Question format is the single-ECG ECG-QA format evaluated in this work. Attribute type is the broad ECG-QA metadata category used to generate the question. n is the number of evaluated question instances in the corresponding question-format and attribute-type subset. EM is exact-match accuracy over the complete predicted answer set. μP , μR , and $\mu F1$ are micro-averaged option-level precision, recall, and F1.

7 MIMIC-IV-ECG ECG-QA performance by question format and attribute

Question format	Attribute type	n	EM	$\mu\mathbf{P}$	$\mu\mathbf{R}$	$\mu\mathbf{F1}$
Single-Verify	SCP code	23,950	81.16	81.16	81.16	81.16
Single-Verify	Heart axis	240	81.67	81.67	81.67	81.67
Single-Verify	Stage of infarction	180	78.33	78.33	78.33	78.33
Single-Verify	Noise	148	79.05	79.05	79.05	79.05
Single-Verify	Numeric feature	1,020	79.22	79.22	79.22	79.22
Single-Choose	SCP code	21,971	68.43	75.41	71.82	73.57
Single-Choose	Heart axis	38	50.00	55.26	52.50	53.85
Single-Choose	Stage of infarction	18	77.78	77.78	77.78	77.78
Single-Choose	Numeric feature	104	60.58	60.58	60.58	60.58
Single-Query	SCP code	54,167	19.25	55.34	48.35	51.61
Single-Query	Heart axis	200	64.50	64.50	64.50	64.50
Single-Query	Stage of infarction	250	44.80	44.80	44.80	44.80
Single-Query	Noise	27	7.41	12.50	12.90	12.70
Single-Query	Numeric feature	7,617	20.90	69.54	66.92	68.21

Table 11: MIMIC-IV-ECG ECG-QA performance stratified by question format and ECG-QA attribute type. The question format is the single-ECG ECG-QA format evaluated in this work. Attribute type is the broad ECG-QA metadata category used to generate the question. n is the number of evaluated question instances in the corresponding question-format and attribute-type subset. EM is exact-match accuracy over the complete predicted answer set. $\mu\mathbf{P}$, $\mu\mathbf{R}$, and $\mu\mathbf{F1}$ are micro-averaged option-level precision, recall, and F1.

8 ECG-QA Single-Verify performance by attribute

Table 12: ECG-QA Single-Verify performance by attribute for PTB-XL and MIMIC-IV-ECG. EM is reported as percent.

Attribute	MIMIC		PTB-XL	
	<i>n</i>	EM	<i>n</i>	EM
SCP-code diagnostic statements				
2:1 AV block	60	93.3	—	—
2:1 sinoatrial block	42	92.9	—	—
3:1 AV block	60	70.0	—	—
4:1 AV block	60	73.3	—	—
aberrant supraventricular complexes	60	75.0	—	—
aberrant ventricular complex	60	90.0	—	—
abnormal qrs	—	—	60	71.7
abnormal R wave progression	60	83.3	—	—
abnormal ventricular conduction pathways	60	73.3	—	—
accelerated idioventricular rhythm	60	88.3	—	—
accelerated junctional rhythm	80	65.0	—	—
any diagnostic symptoms	120	75.8	120	80.8
any form-related symptoms	60	80.0	60	66.7
any kind of abnormal symptoms	60	81.7	60	71.7
any rhythm-related symptoms	60	75.0	60	65.0
atrial arrhythmia	12	100.0	—	—
atrial bigeminy	60	66.7	—	—
atrial couplet	60	65.0	—	—
atrial fibrillation	60	90.0	60	86.7
atrial flutter	80	63.7	36	94.4
atrial premature complex	—	—	60	71.7
atrial tachycardia	60	91.7	—	—
atypical left bundle branch block	60	91.7	—	—
atypical right bundle branch block	60	85.0	—	—
AV dissociation	60	76.7	—	—
AV sequential pacemaker	60	90.0	—	—
bi-atrial overload/enlargement	80	68.8	—	—
bigeminal pattern (unknown origin, supraventricular, or ventricular)	—	—	50	80.0
Biventricular hypertrophy	60	91.7	—	—
Broad R wave in lateral leads	60	86.7	—	—
complete left bundle branch block	60	93.3	60	100.0
complete right bundle branch block	60	98.3	61	96.7
conduction disturbance	60	71.7	60	83.3
Deep S wave	581	91.7	—	—
dextrocardia	56	78.6	—	—
digitalis effect	—	—	67	76.1
dual chamber electronic pacing	60	98.3	—	—
early R wave transition	60	95.0	—	—
early repolarization	7	85.7	—	—
ectopic atrial bradycardia	60	88.3	—	—
ectopic atrial rhythm	60	73.3	—	—
ectopic atrial tachycardia	60	91.7	—	—
electronic atrial pacing	60	83.3	—	—
extreme tachycardia	59	93.2	—	—
first degree AV block	60	83.3	63	66.7
fusion complexes	60	86.7	—	—
high amplitude T-waves	60	91.7	—	—
high grade AV block	58	77.6	—	—
high P-voltages	426	82.4	—	—
high QRS voltage	600	86.8	145	86.2
hypertrophy	60	83.3	60	88.3
idioventricular rhythm	60	86.7	—	—
incomplete left bundle branch block	60	91.7	51	92.2
incomplete right bundle branch block	60	88.3	69	71.0
intermittent second degree AV block	60	80.0	—	—
intraventricular conduction disturbance	60	88.3	—	—
inverted T-waves	60	70.0	621	82.0
ischemic in anterior leads	—	—	32	84.4
ischemic in anterolateral leads	—	—	80	61.3
ischemic in anteroseptal leads	—	—	63	84.1
ischemic in inferior leads	—	—	66	74.2
ischemic in inferolateral leads	—	—	67	62.7
ischemic in lateral leads	—	—	55	65.5
ischemic ST-T changes	60	85.0	—	—
ischemic ST-T changes in anterior leads	80	66.2	—	—
ischemic ST-T changes in anteroseptal leads	80	68.8	—	—
ischemic ST-T changes in diffuse leads	60	65.0	—	—
ischemic ST-T changes in inferior leads	80	62.5	—	—
ischemic ST-T changes in inferoseptal leads	80	77.5	—	—
ischemic ST-T changes in lateral leads	14	71.4	—	—
ischemic ST-T changes in septal leads	80	76.2	—	—
junctional bradycardia	50	80.0	—	—
junctional rhythm	80	67.5	—	—

Attribute	MIMIC		PTB-XL	
	<i>n</i>	EM	<i>n</i>	EM
junctional tachycardia	60	93.3	—	—
late R wave transition	30	86.7	—	—
left anterior fascicular block	80	77.5	68	82.4
left atrial overload/enlargement	80	43.8	80	55.0
left posterior fascicular block	60	85.0	60	98.3
Left ventricular hypertrophy	80	65.0	80	78.8
long QT interval	60	88.3	52	78.8
low amplitude t-wave	—	—	636	78.0
low QRS voltage	1260	93.7	—	—
low qrs voltages in the frontal and horizontal leads	—	—	126	88.1
low R	438	75.3	—	—
Mobitz type 1 second-degree AV block	60	73.3	—	—
Mobitz type 2 second-degree AV block	49	83.7	—	—
multifocal premature ventricular complexes	60	88.3	—	—
Myocardial infarction	60	85.0	60	78.3
Myocardial infarction in anterior leads	80	71.2	80	67.5
Myocardial infarction in anterolateral leads	80	63.7	76	78.9
Myocardial infarction in anteroseptal leads	80	57.5	80	66.2
Myocardial infarction in inferior leads	80	67.5	80	75.0
Myocardial infarction in inferolateral leads	80	73.8	79	68.4
myocardial infarction in inferoposterior leads	—	—	7	71.4
myocardial infarction in inferoposterolateral leads	—	—	7	85.7
Myocardial infarction in lateral leads	80	66.2	49	81.6
Myocardial infarction in posterior leads	80	75.0	7	85.7
Myocardial infarction in septal leads	80	72.5	—	—
non-diagnostic T abnormalities	1260	87.2	743	79.3
non-specific intraventricular conduction disturbance (block)	—	—	63	76.2
non-specific ischemic	—	—	80	70.0
non-specific ST changes	1080	77.6	456	77.4
non-specific ST depression	198	85.9	656	80.5
non-specific st elevation	—	—	24	83.3
non-specific t abnormality	1260	78.4	—	—
non-specific T-wave changes	1080	77.9	684	75.4
non-sustained ventricular tachycardia	60	93.3	—	—
Normal ECG	60	85.0	60	86.7
normal functioning artificial pacemaker	—	—	48	100.0
notched P wave	49	83.7	—	—
P wave abnormality	60	73.3	—	—
pacemaker activity	60	85.0	—	—
pacemaker rhythm	60	96.7	—	—
paired ventricular premature complexes	60	80.0	—	—
paroxysmal idioventricular rhythm	60	73.3	—	—
poor R wave progression	420	76.0	—	—
premature atrial complexes	60	73.3	—	—
premature ventricular complexes	60	85.0	—	—
premature ventricular interpolated complexes	60	83.3	—	—
Prolonged PR interval	60	68.3	60	81.7
prolonged QRS duration	60	85.0	—	—
Q waves present	1039	80.2	647	83.2
QRS changes in anteroseptal leads	60	76.7	—	—
QTc prolongation	60	73.3	—	—
rapid ventricular response	60	88.3	—	—
regular rhythm	60	76.7	—	—
repolarization abnormality	1260	81.0	—	—
Reversed R wave progression	52	82.7	—	—
right atrial overload/enlargement	60	85.0	58	91.4
right ventricular conduction delay	59	71.2	—	—
Right ventricular hypertrophy	80	70.0	—	—
rSr' type in V1 or V2	60	95.0	—	—
S1 S2 S3 type QRS pattern	36	91.7	—	—
second degree AV block	60	78.3	—	—
second-degree SA block type II	60	68.3	—	—
Short PR interval	60	70.0	—	—
short QT interval	60	78.3	—	—
short QTc	30	86.7	—	—
Significant repolarization change	60	73.3	—	—
sinus arrhythmia	60	60.0	60	71.7
sinus bradycardia	60	95.0	60	83.3
sinus pause	60	71.7	—	—
sinus rhythm	60	80.0	60	76.7
sinus tachycardia	80	92.5	60	98.3
slow ventricular response	60	91.7	—	—
ST Depression	1260	84.0	—	—
ST Elevation	1020	77.3	—	—
ST(-T) change	1260	79.1	—	—
st/t change	—	—	60	78.3
subendocardial injury in anterolateral leads	—	—	61	88.5
subendocardial injury in anteroseptal leads	—	—	63	88.9
subendocardial injury in inferolateral leads	—	—	6	83.3
subendocardial injury in lateral leads	—	—	13	84.6
supraventricular bigeminy BIGU bigeminal pattern	60	60.0	—	—

Attribute	MIMIC		PTB-XL	
	<i>n</i>	EM	<i>n</i>	EM
supraventricular extrasystoles	60	61.7	–	–
supraventricular premature complex	60	88.3	–	–
supraventricular rhythm	60	71.7	–	–
supraventricular tachycardia	80	76.2	18	77.8
T-wave abnormality	1260	79.9	150	80.7
tall R wave in V1 or V2	60	93.3	–	–
tall R wave in V5 or V6	57	86.0	–	–
third degree AV block	60	81.7	–	–
TU fusion	52	78.8	–	–
unclassified aberrantly conducted complexes	60	78.3	–	–
uncontrolled ventricular response	60	98.3	–	–
ventricular bigeminy	60	86.7	–	–
ventricular couplet	60	80.0	–	–
ventricular escape rhythm	24	79.2	–	–
ventricular premature complex	–	–	60	88.3
ventricular trigeminy	60	95.0	–	–
ventricular-paced complexes or rhythm	60	95.0	–	–
Voltage criteria (QRS) for left ventricular hypertrophy	60	88.3	180	88.3
wandering pacemaker	36	77.8	–	–
wide QRS tachycardia	60	95.0	–	–
Wolff-Parkinson type A	6	83.3	–	–
Wolff-Parkinson type B	50	82.0	–	–
Wolff-Parkinson-White syndrome	60	85.0	–	–
Heart axis				
extreme axis deviation	60	90.0	57	93.0
left axis deviation	60	66.7	60	80.0
normal heart axis	60	76.7	60	93.3
right axis deviation	60	93.3	60	93.3
Stage of infarction				
early stage of myocardial infarction	60	83.3	59	89.8
middle stage of myocardial infarction	60	80.0	60	80.0
old stage of myocardial infarction	60	71.7	49	77.6
Noise				
any kind of noises	–	–	780	65.9
baseline drift	–	–	780	66.7
baseline wander	148	79.1	–	–
burst noise	–	–	777	68.1
electrodes problems	–	–	66	81.8
static noise	–	–	780	69.1
Numeric features				
above the normal range of p duration	60	68.3	60	68.3
above the normal range of pr interval	60	81.7	60	83.3
above the normal range of qrs duration	60	78.3	60	78.3
above the normal range of qt corrected	60	70.0	60	75.0
above the normal range of qt interval	60	81.7	60	70.0
above the normal range of rr interval	60	91.7	60	98.3
below the normal range of pr interval	60	66.7	60	61.7
below the normal range of qrs duration	60	85.0	60	85.0
below the normal range of qt corrected	60	78.3	60	81.7
below the normal range of qt interval	60	81.7	60	83.3
below the normal range of rr interval	60	100.0	60	93.3
within the normal range of p duration	60	73.3	60	65.0
within the normal range of pr interval	60	68.3	60	71.7
within the normal range of qrs duration	60	78.3	60	63.3
within the normal range of qt corrected	60	70.0	60	73.3
within the normal range of qt interval	60	73.3	60	76.7
within the normal range of rr interval	60	100.0	60	95.0
Extra systole				
any kind of extra systoles	–	–	50	82.0
extrasystoles	–	–	60	81.7
supraventricular extrasystoles	–	–	60	66.7
ventricular extrasystoles	–	–	60	80.0

9 Fields used for structured-record QA generation

Table 13: Dataset-specific fields used to construct structured per-study records for synthetic QA generation.

Field	Dataset	Comment
Patient context		
Patient age	UK Biobank, MIMIC-IV-ECG, EchoNext, PTB-XL	MIMIC-IV-ECG uses approximate anchor age.
Patient sex	UK Biobank, MIMIC-IV-ECG, EchoNext, PTB-XL	–
BMI	UK Biobank	–
Height	UK Biobank	–
Body surface area	UK Biobank	Used for indexed CMR reference ranges.
Chest pain or discomfort	UK Biobank	–
Shortness of breath on level ground	UK Biobank	–
Race or ethnicity	EchoNext	–
ECG measurements		
ECG heart rate	MIMIC-IV-ECG, UK Biobank	In MIMIC-IV-ECG, derived from RR interval if heart rate is unavailable.
Ventricular rate	EchoNext	–
Atrial rate	EchoNext	–
RR interval	MIMIC-IV-ECG, UK Biobank	–
PP interval	MIMIC-IV-ECG, UK Biobank	–
PR interval	MIMIC-IV-ECG, EchoNext	–
PQ interval	MIMIC-IV-ECG, UK Biobank	Derived from P onset and Q onset when needed.
P duration	MIMIC-IV-ECG, UK Biobank	Derived from P onset and P offset when needed.
QRS duration	MIMIC-IV-ECG, UK Biobank, EchoNext	In MIMIC-IV-ECG, derived from QRS onset and QRS end when needed.
QT interval	MIMIC-IV-ECG, UK Biobank	Used to derive QTc when needed.
QTc interval	MIMIC-IV-ECG, UK Biobank, EchoNext	In MIMIC-IV-ECG and UK Biobank, derived from QT and RR when needed.
P axis	MIMIC-IV-ECG, UK Biobank	–
R axis	MIMIC-IV-ECG, UK Biobank, PTB-XL	In PTB-XL, heart-axis labels are mapped to unified R-axis categories. MID: normal; LAD and ALAD: leftward; RAD and ARAD: rightward.
T axis	MIMIC-IV-ECG, UK Biobank	–
ECG report-derived fields		
ECG diagnosis and machine-report diagnosis text	MIMIC-IV-ECG, UK Biobank	ECG diagnosis is retained as text. For MIMIC-IV-ECG, diagnosis text is built from machine-measurement report columns.
German ECG report	PTB-XL	Original PTB-XL report text.
English ECG report	PTB-XL	English translation of the PTB-XL report generated using GPT-OSS-20B.
SCP-ECG diagnostic codes	PTB-XL	SCP codes are decoded into textual diagnostic descriptions.
Infarction stage	PTB-XL	–
CMR-derived fields		
LV end diastolic volume	UK Biobank	–
LV end systolic volume	UK Biobank	–
LV stroke volume	UK Biobank	–
LV myocardial mass	UK Biobank	–
LV ejection fraction	UK Biobank	–
LV cardiac output	UK Biobank	–
Cardiac index	UK Biobank	–
LV mean myocardial wall thickness global	UK Biobank	–
LV circumferential strain global	UK Biobank	–
LV longitudinal strain global	UK Biobank	–
LV radial strain global	UK Biobank	–
RV end diastolic volume	UK Biobank	–

Field	Dataset	Comment
RV end systolic volume	UK Biobank	–
RV stroke volume	UK Biobank	–
RV ejection fraction	UK Biobank	–
LA maximum volume	UK Biobank	–
LA minimum volume	UK Biobank	–
LA stroke volume	UK Biobank	–
LA ejection fraction	UK Biobank	–
RA maximum volume	UK Biobank	–
RA minimum volume	UK Biobank	–
RA stroke volume	UK Biobank	–
RA ejection fraction	UK Biobank	–
Ascending aorta distensibility	UK Biobank	–
Descending aorta distensibility	UK Biobank	–
EchoNext echocardiographic categorical fields		
Aortic stenosis	EchoNext	Severity category.
Aortic regurgitation	EchoNext	Severity category.
Mitral regurgitation	EchoNext	Severity category.
Tricuspid regurgitation	EchoNext	Severity category.
Pulmonary regurgitation	EchoNext	Severity category.
Right ventricular systolic function	EchoNext	Qualitative systolic function category.
Pericardial effusion	EchoNext	Presence and size category.
EchoNext echocardiographic numeric fields		
Interventricular septal thickness	EchoNext	–
Left ventricular posterior wall thickness	EchoNext	–
Left ventricular systolic function	EchoNext	Derived from LVEF value in percent.
EchoNext binary structural heart disease fields		
Left ventricular ejection fraction is 45% or lower	EchoNext	Binary descriptor derived from LVEF.
Left ventricular wall thickness is at least 1.3 cm	EchoNext	Binary descriptor derived from LV wall thickness.
Aortic stenosis is moderate or greater	EchoNext	Binary descriptor derived from aortic stenosis severity.
Aortic regurgitation is moderate or greater	EchoNext	Binary descriptor derived from aortic regurgitation severity.
Mitral regurgitation is moderate or greater	EchoNext	Binary descriptor derived from mitral regurgitation severity.
Tricuspid regurgitation is moderate or greater	EchoNext	Binary descriptor derived from tricuspid regurgitation severity.
Pulmonary regurgitation is moderate or greater	EchoNext	Binary descriptor derived from pulmonary regurgitation severity.
Right ventricular systolic dysfunction is moderate or greater	EchoNext	Binary descriptor derived from right ventricular systolic function.
Pericardial effusion is moderate or large	EchoNext	Binary descriptor derived from pericardial effusion size.
Pulmonary artery systolic pressure is at least 45 mmHg	EchoNext	Binary descriptor derived from PASP.
Tricuspid regurgitation max velocity is at least 3.2 m/s	EchoNext	Binary descriptor derived from maximum TR jet velocity.

10 Clinical categorization rules

Table 14: Clinical categorisation rules used to convert numeric measurements into textual descriptors before synthetic QA generation.

Measurement	Categorisation rule
ECG measurements	
ECG heart rate [bpm], ventricular rate, atrial rate	≤ 50 : marked bradycardia; > 50 to ≤ 60 : bradycardia; > 60 to ≤ 100 : normal; > 100 to ≤ 120 : mild tachycardia; > 120 : marked tachycardia.
RR interval [ms], PP interval [ms]	Convert interval to rate as $60000/\text{interval}$. Rate > 120 : markedly short; > 100 to ≤ 120 : short; ≥ 60 to ≤ 100 : normal; ≥ 50 to < 60 : prolonged; < 50 : markedly prolonged.
P duration [ms]	If missing, derive as P offset – P onset. Then categorise as: < 120 : normal; ≥ 120 : prolonged.
PQ interval [ms], PR interval [ms]	If PQ is missing, derive as Q onset – P onset. Then categorise as: < 120 : short; 120–200: normal; > 200 : prolonged.
QRS duration [ms]	< 110 : normal; ≥ 110 to < 120 : mildly prolonged; ≥ 120 : prolonged.
QTC interval [ms]	If QTC is unavailable, derive it from QT and RR using Bazett correction: $QTc = QT/\sqrt{RR/1000}$. Female: ≤ 470 : normal; > 470 to ≤ 490 : borderline; > 490 : prolonged. Male or unavailable sex: ≤ 450 : normal; > 450 to ≤ 480 : borderline; > 480 : prolonged.
P axis [degrees]	< 0 : leftward; 0–75: normal; > 75 : rightward.
R axis [degrees]	< -30 : leftward; -30 –90: normal; > 90 : rightward.
T axis [degrees]	< -15 : leftward; -15 to < 15 : borderline; 15–75: normal; > 75 to ≤ 105 : borderline; > 105 : rightward.
Ventricular function and cardiac output	
LV ejection fraction, left ventricular systolic function	Female: < 30 : severely reduced; ≥ 30 to < 40 : moderately reduced; ≥ 40 to < 52 : mildly reduced; 52–79: normal; > 79 : hyperdynamic. Male or unavailable sex: < 30 : severely reduced; ≥ 30 to < 40 : moderately reduced; ≥ 40 to < 49 : mildly reduced; 49–79: normal; > 79 : hyperdynamic.
RV ejection fraction	Female: < 30 : severely reduced; ≥ 30 to < 40 : moderately reduced; ≥ 40 to < 46 : mildly reduced; 46–74: normal; > 74 : hyperdynamic. Male or unavailable sex: < 30 : severely reduced; ≥ 30 to < 40 : moderately reduced; ≥ 40 to < 42 : mildly reduced; 42–72: normal; > 72 : hyperdynamic.
cardiac index [litres/min/m ²]	< 2.2 : low output; 2.2 to ≤ 2.5 : borderline low; > 2.5 to ≤ 4.0 : normal; > 4.0 to ≤ 4.5 : borderline high; > 4.5 : high output.
LV cardiac output	If body surface area is available, first compute cardiac index and apply the cardiac-index rule. Otherwise, use absolute ranges. Male 3.4–7.8: normal; female 2.7–6.3: normal; unknown sex 2.7–7.8: normal. Below range: low output; above range: high output.
LV mean myocardial wall thickness global	< 5.0 : thinned; 5.0–11.0: normal; > 11.0 : hypertrophied.
Echocardiographic LV wall thickness [cm], interventricular septal thickness, left ventricular posterior wall thickness	< 0.6 : thinned. Male: ≤ 1.0 : normal; > 1.0 to ≤ 1.3 : mildly increased; > 1.3 to ≤ 1.6 : moderately increased; > 1.6 : severely increased. Female or unknown sex: ≤ 0.9 : normal; > 0.9 to ≤ 1.2 : mildly increased; > 1.2 to ≤ 1.5 : moderately increased; > 1.5 : severely increased.
CMR ventricular volumes, stroke volume, and mass	
LV end diastolic volume	If body surface area is available, use indexed range: male 50–108: normal; female 50–96: normal. Otherwise use absolute range: male 95–215: normal; female 78–167: normal. Below range: small; above range: dilated.
LV end systolic volume	If body surface area is available, use indexed range: male 11–47: normal; female 10–40: normal. Otherwise use absolute range: male 23–100: normal; female 18–73: normal. Below range: small; above range: dilated.
LV stroke volume	If body surface area is available, use indexed range: male 33–72: normal; female 33–64: normal. Otherwise use absolute range: male 60–135: normal; female 47–99: normal. Below range: reduced; above range: increased.
LV myocardial mass	If body surface area is available, use indexed range: male 39–85: normal; female 30–68: normal. Otherwise use absolute range: male 73–171: normal; female 61–121: normal. Below range: decreased; above range: increased.
RV end diastolic volume	If body surface area is available, use indexed range: male 53–123: normal; female 48–104: normal. Otherwise use absolute range: male 105–258: normal; female 84–193: normal. Below range: small; above range: dilated.
RV end systolic volume	If body surface area is available, use indexed range: male 17–59: normal; female 13–48: normal. Otherwise use absolute range: male 34–123: normal; female 22–86: normal. Below range: small; above range: dilated.
RV stroke volume	If body surface area is available, use indexed range: male 28–75: normal; female 29–66: normal. Otherwise use absolute range: male 57–141: normal; female 44–116: normal. Below range: reduced; above range: increased.
CMR atrial volumes and function	

Measurement	Categorisation rule
LA maximum volume	If body surface area is available, use indexed range: male 17–59: normal; female 17–61: normal. Otherwise use absolute range: male 31–112: normal; female 28–100: normal. Below range: small; above range: dilated.
LA minimum volume	If body surface area is available, use indexed range: male 3–24: normal; female 4–23: normal. Otherwise use absolute range: male 6–44: normal; female 7–38: normal. Below range: small; above range: dilated.
LA stroke volume	If body surface area is available, use indexed range: male 10–34: normal; female 10–34: normal. Otherwise use absolute range: male 21–67: normal; female 21–62: normal. Below range: reduced; above range: increased.
LA ejection fraction	Male 43–75: normal; female 47–75: normal. Below range: reduced; above range: hyperdynamic.
RA maximum volume	If body surface area is available, use indexed range: male 32–79: normal; female 31–69: normal. Otherwise use absolute range: male 59–158: normal; female 49–122: normal. Below range: small; above range: dilated.
RA minimum volume	If body surface area is available, use indexed range: male 10–42: normal; female 8–32: normal. Otherwise use absolute range: male 16–84: normal; female 11–55: normal. Below range: small; above range: dilated.
RA stroke volume	If body surface area is available, use indexed range: male 14–46: normal; female 14–42: normal. Otherwise use absolute range: male 26–90: normal; female 23–71: normal. Below range: reduced; above range: increased.
RA ejection fraction	Male 34–74: normal; female 41–77: normal. Below range: reduced; above range: hyperdynamic.
Global strain	
LV circumferential strain global	Male –27.2 to –14.6: normal; female –29.2 to –16.2: normal. Value above the upper limit: reduced magnitude; value below the lower limit: increased magnitude.
LV longitudinal strain global	Male –26.1 to –12.7: normal; female –28.7 to –14.2: normal. Value above the upper limit: reduced magnitude; value below the lower limit: increased magnitude.
LV radial strain global	Male and female 20.0–55.0: normal. Below range: reduced magnitude; above range: increased magnitude.
Aortic distensibility	
Ascending aorta distensibility	Male age ≤ 54 : 0.7–5.1: normal; male age 55–64: 0.0–4.2: normal; male age > 64 : 0.0–2.4: normal. Female age ≤ 54 : 0.5–5.7: normal; female age 55–64: 0.0–3.9: normal; female age > 64 : 0.0–2.7: normal. Below range: decreased; above range: increased.
Descending aorta distensibility	Male age ≤ 54 : 1.6–6.0: normal; male age 55–64: 0.7–5.1: normal; male age > 64 : 0.4–4.0: normal. Female age ≤ 54 : 1.6–6.0: normal; female age 55–64: 0.6–4.6: normal; female age > 64 : 0.4–3.6: normal. Below range: decreased; above range: increased.
EchoNext binary echocardiography rules	
$LVEF \leq 45\%$	Binary descriptor: yes if left ventricular ejection fraction is 45% or lower; otherwise no.
$LV \text{ wall thickness} \geq 1.3 \text{ cm}$	Binary descriptor: yes if left ventricular wall thickness is at least 1.3 cm; otherwise no.
$Aortic \text{ stenosis} \geq moderate$	Binary descriptor: yes if aortic stenosis is moderate or greater; otherwise no.
$Aortic \text{ regurgitation} \geq moderate$	Binary descriptor: yes if aortic regurgitation is moderate or greater; otherwise no.
$Mitral \text{ regurgitation} \geq moderate$	Binary descriptor: yes if mitral regurgitation is moderate or greater; otherwise no.
$Tricuspid \text{ regurgitation} \geq moderate$	Binary descriptor: yes if tricuspid regurgitation is moderate or greater; otherwise no.
$Pulmonary \text{ regurgitation} \geq moderate$	Binary descriptor: yes if pulmonary regurgitation is moderate or greater; otherwise no.
$RV \text{ systolic dysfunction} \geq moderate$	Binary descriptor: yes if right ventricular systolic dysfunction is moderate or greater; otherwise no.
$Pericardial \text{ effusion moderate or large}$	Binary descriptor: yes if pericardial effusion is moderate or large; otherwise no.
$PASP \geq 45 \text{ mmHg}$	Binary descriptor: yes if pulmonary artery systolic pressure is at least 45 mmHg; otherwise no.
$TR \text{ maximum velocity} \geq 3.2 \text{ m/s}$	Binary descriptor: yes if tricuspid regurgitation maximum velocity is at least 3.2 m/s; otherwise no.

11 ECG-LLM Chat Template

During ECG-LLM instruction tuning and evaluation, each example was serialized as a chat conversation with a fixed system message, a user message containing the ECG token placeholder and clinical question, and an assistant message containing the target answer. The ECG token placeholder denotes the projected ECG-token sequence inserted into the LLM input embedding stream.

ECG-LLM chat template

System:

You are a cardiology assistant. The user will provide ECG data after the phrase "Here is the
↪ ECG:" followed by a clinical question. Answer the question directly using only findings
↪ supported by the ECG and available patient context. Do not invent unsupported findings,
↪ reason only on the provided ECG.

User:

Here is the ECG: <ecg_tokens>
<question_tokens>

Assistant:

<answer_tokens>

12 Question Answer Generation Prompts

The following prompts were used to generate question-answer pairs for training. Dataset-specific prompt variants were used for UK Biobank, EchoNext, MIMIC-IV-ECG, and PTB-XL according to the structured fields available in each cohort. For UK Biobank and EchoNext, we additionally generated reasoning-style QA pairs. In these examples, the generator returned separate `reasoning` and `answer` fields. During training, the model input was the generated question, and the target response was formed by concatenating the reasoning field with the final answer field.

Prompt for standard QA generation with UK Biobank data

Create N question-answer pairs from this text for LLM training in English.

BEHAVIOR:

- Simulate an interaction via QA pairs between a General Practitioner (GP) and a perfect AI system for ECG interpretation.
- The GP asks questions about this ECG, and the AI answers.
- The GP sees only the ECG and asks clinically meaningful, open-ended questions about this patient to obtain a comprehensive understanding of the patient's cardiac condition and abnormalities.
- The AI may internally use all provided information from the text.
- The AI must always answer as if the insights come purely from ECG interpretation.
- Answers should use as much provided information as possible, combining several fields to provide a comprehensive picture of the patient's cardiac status and abnormalities.
- Do not invent diagnoses, numbers, thresholds, or advice not supported by the provided fields.
- These QA pairs are intended for LLM fine-tuning. Across the full set of generated pairs, they should collectively reflect the patient's overall cardiac state, major abnormalities, and all clinically relevant supported findings in the text.

QUESTIONS (GP):

- Questions must sound like a GP asking about the ECG of this specific patient.
- Questions must not state or assume specific ECG diagnoses, measurements, or values.
- Do not ask about information that is not present in the text.
- Do not include the findings in the question. For example, do not ask: "If sinus bradycardia is present...", "Given LV dilatation...", or "If EF is reduced...".
- Instead, ask the AI to interpret, describe, and highlight abnormalities and to identify all relevant findings.
- Questions may reflect uncertainty or clinical curiosity, for example: "Does the ECG suggest ...?", "Could this ECG indicate ...?", "Do you see any signs of ...?", or "Is there evidence of ...?".

ANSWERS (AI):

- The AI has access to all information in the text, but must formulate answers as if the insights come purely from ECG interpretation.
- Answers must contain several clear sentences, with a direct and unambiguous answer to the question.
- Answers must include as much relevant information as possible about the patient's cardiac condition and abnormalities, combining several fields where appropriate.
- All answers must be directly supported by the provided fields and must never contradict any given measurements.
- Never mention anything not present in the provided text.

OUTPUT FORMAT (STRICT):

Return only valid JSON in this exact format:

```
[
  {
    "question": "Question 1?",
    "answer": "Answer 1."
  },
]
```

```

    {{
      "question": "Question 2?",
      "answer": "Answer 2."
    }}
  ]

```

```

Text:
{text}

```

Prompt for reasoning-style QA generation with UK Biobank data

Create N complex reasoning-style question-answer examples from this text for LLM training in
 ↪ English.

BEHAVIOR:

- Simulate an interaction via QA pairs between a General Practitioner (GP) and a perfect AI
 ↪ system for ECG interpretation.
- The GP asks questions about this ECG, and the AI answers.
- The GP sees only the ECG and asks clinically meaningful, open-ended questions about this
 ↪ patient to obtain a comprehensive understanding of the patient's cardiac condition and
 ↪ abnormalities.
- The AI may internally use all provided information from the text.
- The AI must always answer as if the insights come purely from ECG interpretation.
- Answers should use as much provided information as possible, combining several fields to
 ↪ provide a comprehensive picture of the patient's cardiac status and abnormalities.
- Do not invent diagnoses, numbers, thresholds, or advice not supported by the provided
 ↪ fields.
- These examples are intended for LLM fine-tuning. Across the full set of generated examples,
 ↪ they should collectively reflect the patient's overall cardiac state, major
 ↪ abnormalities, and all clinically relevant supported findings in the text.

QUESTIONS (GP):

- Questions must sound like a GP asking about the ECG of this specific patient.
- Questions must be clinically meaningful and require multi-step interpretation.
- Questions must not state or assume specific ECG diagnoses, measurements, or values.
- Do not ask about information that is not present in the text.
- Do not include the findings in the question. For example, do not ask: "If sinus bradycardia
 ↪ is present...", "Given LV dilatation...", or "If EF is reduced...".
- Instead, ask the AI to interpret, describe, and highlight abnormalities and to identify all
 ↪ relevant findings.
- Questions may reflect uncertainty or clinical curiosity, for example: "Does the ECG suggest
 ↪ ...?", "Could this ECG indicate ...?", "Do you see any signs of ...?", or "Is there
 ↪ evidence of ...?".

ANSWERS (AI):

- The AI has access to all information in the text, but must formulate answers as if the
 ↪ insights come purely from ECG interpretation.
- Each example must contain:
 1. a challenging ECG-focused clinical question,
 2. a concise reasoning field that links the available findings step by step,
 3. a final answer that directly addresses the question.
- The reasoning must be clinically grounded and directly supported by the provided fields.
- The final answer must be concise, direct, and unambiguous.
- All reasoning and answers must be directly supported by the provided fields and must never
 ↪ contradict any given measurements.
- Never mention anything not present in the provided text.

OUTPUT FORMAT (STRICT):

Return only valid JSON in this exact format:

```

[
  {{
    "question": "Complex question about this ECG?",

```

```

    "reasoning": "Step 1: ... Step 2: ... Step 3: ...",
    "answer": "Final answer based on the reasoning."
  }},
  {{
    "question": "Another complex question about this ECG?",
    "reasoning": "Step 1: ... Step 2: ... Step 3: ...",
    "answer": "Final answer based on the reasoning."
  }}
]

Text:
{text}

```

Prompt for QA generation with MIMIC-IV-ECG data

Create N question-answer pairs from this text for LLM training in English.

BEHAVIOR:

- Simulate an interaction between a General Practitioner (GP) and a PERFECT AI for ECG
 - ↪ interpretation.
- The GP asks questions about this ECG, and the PERFECT AI answers.
- The GP sees only the ECG and asks clinically meaningful, open-ended questions about THIS
 - ↪ patient to obtain a comprehensive understanding of the patient's heart condition and
 - ↪ abnormalities.
- PERFECT AI for ECG interpretation may internally use and reason with ALL provided
 - ↪ information from the text (ECG phenotypes: + Patient demographics),
 - but MUST answer as if insights come purely from perfect ECG interpretation.
- Questions and answers must be phrased as if inferred from ECG alone. Never mention records,
 - ↪ reports, metadata, or models.
- Do not invent diagnoses, numbers, thresholds, or advice not supported by the provided
 - ↪ fields.
- These QA pairs are intended for LLM fine-tuning. Across the full set of generated pairs,
 - ↪ they should collectively reflect the patient's overall cardiac state, major
 - ↪ abnormalities, and all clinically relevant supported findings in the text.

QUESTIONS (GP):

- Must sound like a GP interpreting and asking questions about the ECG of this specific
 - ↪ patient.
- MUST NOT state or assume specific ECG diagnoses, measurements, or values in the question.
- Do not ask about information that is not present in the text.
- Do not include the findings in the question. Instead, ask the AI to interpret, describe,
 - ↪ and highlight abnormalities and to explain all relevant findings.
- May reflect uncertainty or clinical curiosity: "Does the ECG suggest ...?", "Could this ECG indicate...?", "Do you see any signs of...", "Is there evidence of...", etc.

ANSWERS (AI):

- This is a PERFECT AI for ECG interpretation. It has access to all information in the text,
 - ↪ but must formulate answers as if insights come purely from ECG interpretation.
- Answers must be several clear sentences, with a direct and unambiguous answer to the
 - ↪ question.
- Answers must include as much relevant information as possible about the patient's heart
 - ↪ condition and abnormalities, combining several fields.
- All answers must be directly supported by the provided fields and must never contradict any
 - ↪ given measurements.
- Never mention anything not present in the provided text.

OUTPUT FORMAT (STRICT):

Return ONLY valid JSON in this exact format:

```

[
  {{
    "question": "Question 1?",
    "answer": "Answer 1."
  }}
]

```

```

    }},
    {{
      "question": "Question 2?",
      "answer": "Answer 2."
    }}
  ]

```

```

Text:
{text}

```

Prompt for QA generation with PTB-XL data

Create N question-answer pairs from this text for LLM training in English.

BEHAVIOR:

- Simulate an interaction between a General Practitioner (GP) and a PERFECT AI for ECG
 - ↪ interpretation.
- The GP asks questions about this ECG, and the PERFECT AI answers.
- The GP sees only the ECG and asks clinically meaningful, open-ended questions about THIS patient to obtain a comprehensive understanding of the patient's heart condition and
 - ↪ abnormalities.
- PERFECT AI for ECG interpretation may internally use and reason with ALL provided
 - ↪ information from the text.
- PERFECT AI for ECG interpretation ALWAYS answer as if insights come purely from ECG
 - ↪ analysis.
- The output QA pairs MUST be ONLY in English.
- Do not invent diagnoses, numbers, thresholds, or advice not supported by the provided
 - ↪ fields.
- These QA pairs are intended for LLM fine-tuning. Across the full set of generated pairs,
 - ↪ they should collectively reflect the patient's overall cardiac state, major
 - ↪ abnormalities, and all clinically relevant supported findings in the text.
- QA pairs must cover all findings from the text.

QUESTIONS (GP):

- Must sound like a GP interpreting and asking questions about the ECG of this specific
 - ↪ patient.
- NEVER include the findings from the text in the question. For example, do NOT say "If sinus
 - ↪ bradycardia is present...", "Given LV dilatation...", or "If EF is reduced...". Instead,
 - ↪ ask the AI to interpret, describe, and highlight abnormalities and to explain all
 - ↪ relevant findings.
- May reflect uncertainty or clinical curiosity: "Does the ECG suggest ...?", "Could this ECG indicate...?", "Do you see any signs of...", "Is there evidence of...", etc.
- Do not ask about information that is not present in the text.

ANSWERS (AI):

- This is a PERFECT AI for ECG interpretation. It must formulate answers as if insights come
 - ↪ purely from ECG interpretation.
- Answers must be several clear sentences, with a direct and unambiguous answer to the
 - ↪ question.
- Answers must include as much relevant information as possible about the patient's heart
 - ↪ condition and abnormalities, combining several fields.
- All answers must be directly supported by the provided fields and must never contradict any
 - ↪ given measurements.
- Never mention anything not present in the provided text.

OUTPUT FORMAT (STRICT):

Return ONLY valid JSON in this exact format:

```

[
  {{
    "question": "Question 1?",
    "answer": "Answer 1."
  }},

```

```

    {{
      "question": "Question 2?",
      "answer": "Answer 2."
    }}
  ]

```

Text:
{text}

Prompt for QA generation with EchoNext data

Create N question-answer pairs from this text for LLM training in English.

BEHAVIOR:

- Simulate an interaction via QA pairs between a General Practitioner (GP) and a PERFECT AI
 - ↪ for ECG interpretation.
- The GP asks questions about this ECG, and the PERFECT AI answers.
- The GP sees only the ECG and asks clinically meaningful, open-ended questions about THIS
 - ↪ patient to obtain a comprehensive understanding of the patient's heart condition and
 - ↪ abnormalities.
- PERFECT AI for ECG interpretation may internally use and reason with ALL provided
 - ↪ information from the text.
- PERFECT AI for ECG interpretation ALWAYS answer as if insights come purely from ECG.
- Answers should cite as much provided information as possible, combining several fields to
 - ↪ provide a comprehensive picture of the patient's cardiac status and abnormalities.
- Do not invent diagnoses, numbers, thresholds, or advice not supported by the provided
 - ↪ fields.
- These QA pairs are intended for LLM fine-tuning. Across the full set of generated pairs,
 - ↪ they should collectively reflect the patient's overall cardiac state, major
 - ↪ abnormalities, and all clinically relevant supported findings in the text.

QUESTIONS (GP):

- Must sound like a GP asking questions about the ECG of this specific patient. GP asks
 - ↪ questions.
- MUST NOT state or assume specific ECG diagnoses, measurements, or values in the question.
- Do not ask about information that is not present in the text.
- Do not include the findings in the question. For example, do NOT say "If sinus bradycardia
 - ↪ is present...", "Given LV dilatation...", or "If EF is reduced...".
- Instead, ask the AI to interpret, describe, and highlight abnormalities and to find all
 - ↪ relevant abnormalities.
- Questions may reflect uncertainty or clinical curiosity: "Does the ECG suggest ...?",
 - ↪ "Could this ECG
 - ↪ indicate...?", "Do you see any signs of...", "Is there evidence of...", etc.

ANSWERS (AI):

- This is a PERFECT AI for ECG interpretation. It has access to all information in the text,
 - ↪ and it must formulate answers as if insights come purely from ECG interpretation.
- Answers must be several clear sentences, with a direct and unambiguous answer to the
 - ↪ question.
- Answers must include as much relevant information as possible about the patient's heart
 - ↪ condition and abnormalities, combining several fields.
- All answers must be directly supported by the provided fields and must never contradict any
 - ↪ given measurements.
- Never mention anything not present in the provided text.

OUTPUT FORMAT (STRICT):

Return ONLY valid JSON in this exact format:

```

[
  {{
    "question": "Question 1?",
    "answer": "Answer 1."
  }},

```

```

    {{
      "question": "Question 2?",
      "answer": "Answer 2."
    }}
  ]

```

Text:
{text}

Prompt for reasoning-style QA generation with EchoNext data

Create N complex reasoning examples from this text that demonstrate chain-of-thought
 ↳ thinking.

BEHAVIOR:

- Simulate an interaction via QA pairs between a General Practitioner (GP) and a PERFECT AI
 ↳ for ECG interpretation.
- The GP asks questions about this ECG, and the PERFECT AI answers.
- The GP sees only the ECG and asks clinically meaningful, open-ended questions about THIS
 ↳ patient to obtain a comprehensive understanding of the patient's heart condition and
 ↳ abnormalities.
- PERFECT AI for ECG interpretation may internally use and reason with ALL provided
 ↳ information from the text.
- PERFECT AI for ECG interpretation ALWAYS answer as if insights come purely from ECG.
- Answers should cite as much provided information as possible, combining several fields to
 ↳ provide a comprehensive picture of the patient's cardiac status and abnormalities.
- Do not invent diagnoses, numbers, thresholds, or advice not supported by the provided
 ↳ fields.
- These QA pairs are intended for LLM fine-tuning. Across the full set of generated pairs,
 ↳ they should collectively reflect the patient's overall cardiac state, major
 ↳ abnormalities, and all clinically relevant supported findings in the text.

QUESTIONS (GP):

- Must sound like a GP asking questions about the ECG of this specific patient. GP asks
 ↳ questions.
- MUST NOT state or assume specific ECG diagnoses, measurements, or values in the question.
- Do not ask about information that is not present in the text.
- Do not include the findings in the question. For example, do NOT say "If sinus bradycardia
 ↳ is present...", "Given LV dilatation...", or "If EF is reduced...".
- Instead, ask the AI to interpret, describe, and highlight abnormalities and to find all
 ↳ relevant abnormalities.
- Questions may reflect uncertainty or clinical curiosity: "Does the ECG suggest ...?",
 ↳ "Could this ECG
 indicate...?", "Do you see any signs of...", "Is there evidence of...", etc.

ANSWERS (AI):

- This is a PERFECT AI for ECG interpretation. It has access to all information in the text,
 ↳ and it must formulate answers as if insights come purely from ECG interpretation.
- Answers must be several clear sentences, with a direct and unambiguous answer to the
 ↳ question.
- Answers must include as much relevant information as possible about the patient's heart
 ↳ condition and abnormalities, combining several fields.
- All answers must be directly supported by the provided fields and must never contradict any
 ↳ given measurements.
- Never mention anything not present in the provided text.

Each example should have:

1. A challenging question that requires step-by-step reasoning
2. Detailed reasoning steps that break down the problem
3. A concise final answer

Return JSON format only:

```
[
  {{
    "question": "Complex question about the text?",
    "reasoning": "Step 1: First, I need to consider...\nStep 2: Then, I analyze...\nStep 3:
    ↪ Finally, I can conclude...",
    "answer": "Final answer based on the reasoning."
  }},
  {{
    "question": "Another complex question?",
    "reasoning": "First, I'll analyze... Next, I need to determine... Based on this
    ↪ analysis...",
    "answer": "Final answer drawn from the reasoning."
  }}
]
```

```
Text:
{text}
```

13 Implementation Details

Parameter	Value
Input	12 leads $\times$ 5,000 samples
Patch size	1×100
ECG tokens	300 before masking
Masking ratio	0.75
Loss	MSE over masked patches
Encoder	ViT, 12 layers, hidden size 768, 12 heads
Decoder	Transformer, 2 layers, hidden size 512, 8 heads
Activation	GELU
Positional embedding	Learnable
Drop-path rate	0.1
Optimiser	AdamW
Learning rate	2×10^{-4}
Weight decay	0.05
Batch size	768
Precision	bfloat16 mixed precision
Training duration	120 epochs
Warm-up	1,000 steps
Minimum LR factor	0.01
Augmentations	Random crop, Fourier surrogate perturbation, jitter, rescaling

Table 15: ECG ViT-MAE pretraining configuration.

Parameter	Value
ECG encoder	Pretrained ViT-MAE, frozen
Language model	Llama-3.1-8B-Instruct, frozen backbone
ECG crop	Random crop to 2,500 time-steps
ECG augmentations	Disabled
Projector	Linear($d_{\text{ecg}}, d_{\text{llm}}$), SiLU, Linear($d_{\text{llm}}, d_{\text{llm}}$)
Aggregation	Identity; all ECG tokens preserved
Projector bias	Disabled
LoRA modules	q.proj , k.proj , v.proj , o.proj
LoRA rank / alpha	128 / 128
LoRA dropout	0.1
Curriculum	First 1,500 steps update only the projector, then projector + LoRA training
Optimiser	Paged AdamW 8-bit
Learning rates	Projector 10^{-4} ; LoRA 10^{-5}
Weight decay	10^{-3}
Schedule	Linear warm-up from $0.1 \times$ to peak LR over 3% of optimiser steps, then cosine decay to $0.3 \times$ peak LR
Precision	bfloat16 mixed precision
Training duration	1 epoch
Batching	4 GPUs, batch size 8/GPU, gradient accumulation 4
Effective batch size	128
Max sequence length	1200 tokens

Table 16: ECG-LLM instruction-tuning configuration.

14 Diagnostic report generation prompts

Model	Diagnostic report generation prompt
ECG-LLM	Review the ECG signal and produce a detailed report on your diagnostic observations, ending with the final diagnosis.
	Describe this ECG in detail and tell if there are any abnormalities present.
	Review this ECG and provide an ECG diagnosis.
	What is wrong with this ECG?
	Provide an ECG report for this tracing.
	What are the key diagnostic findings in this ECG?
	Analyze this ECG carefully and provide a detailed ECG report.
ECG-Chat PULSE-7B MedGemma-4B	What is the summary for this ECG?
	Could you please help me explain my ECG?
	Please write a clinical report based on this ECG image.
	Please interpret this 12-lead ECG image. Provide a concise diagnostic report, including rhythm, rate if assessable, conduction or axis abnormalities if present, ST-T or infarction-related findings if present, and an overall ECG impression. If a finding is uncertain or not clearly visible, say so rather than inventing it.

Table 17: Prompts used for ECG diagnostic summary generation. For ECG-LLM, one question was randomly selected from the ECG-LLM prompt pool for each ECG. For baseline models, the model-specific prompt shown in the table was used.

15 LLM-as-judge prompt for ECG diagnosis evaluation

The following prompts were used for LLM-based evaluation of generated ECG diagnostic answers against ground-truth diagnostic statements.

System prompt for LLM-as-judge evaluation

You are a strict but clinically fair ECG diagnosis evaluator.

Compare MODEL ANSWER against GROUND TRUTH. Treat clinically equivalent wording as correct. Do not
↪ require identical phrasing.
Score correctness, completeness, and specificity independently.

Clinical equivalence rules:

- Borderline ECG, abnormal ECG, nonspecific, possible, probable, normal variant, and age-undetermined
↪ findings are lower-confidence/global findings.
- Do not strongly penalize missing uncertain findings if the main ECG interpretation is correct.
- Penalize missing GT findings under completeness and unsupported extra findings under specificity. Do
↪ not assign 0 if at least one meaningful GT finding is correctly identified. But don't assign 2 if
↪ important GT findings are missing or if the model adds unsupported major findings.

For PTB-XL ground truth: #Added only for PTB-XL evaluation

- Use the Original report as the primary reference.
- Use the English report only as a translation aid.
- Do not treat the English report as an additional independent list of required findings.

Return only valid JSON.

User prompt template for LLM-as-judge evaluation

GROUND TRUTH:
{gt}

MODEL ANSWER:
{pred}

Score independently.

Correctness:
2 = main ECG interpretation is clinically correct or nearly correct.
1 = partially correct: some main findings are correct, but there is an important wrong or missing
↪ conclusion.
0 = mostly incorrect, or the main diagnosis conclusion is wrong.

Completeness:
2 = all or nearly all important GT diagnoses are present.
1 = at least one important GT diagnosis is present, but relevant GT findings are missing.
0 = no clinically meaningful GT diagnosis is present.

Specificity:
2 = nearly all model diagnoses are supported by or compatible with GT.
1 = mixed answer: some supported diagnoses plus unsupported or wrong extras.
0 = mostly unsupported or dominated by wrong diagnoses.

Return JSON exactly:

```
{{
  "correctness": int,
  "completeness": int,
  "specificity": int,
  "total_score": int,
  "notes": "reason"
}}
```