

SAMRI-3D: Adapting SAM2 for 3D MRI Segmentation with Global Volume Tokens

Zhao Wang^{1,*}, Wei Dai¹, Hongfu Sun², Craig Engstrom³, and Shekhar S. Chandra¹

¹School of Electrical Engineering and Computer Science, The University of Queensland, Brisbane, Australia

²School of Engineering, College of Engineering, Science and Environment, University of Newcastle, Australia

³School of Human Movement and Nutrition Sciences, The University of Queensland, Australia

*Corresponding author: zhao.wang1@uq.edu.au

Abstract

Foundation models such as Segment Anything Model 2 (SAM2) have transformed natural-image and video segmentation, and recent work has begun adapting them to medical imaging. These adaptations, however, are largely general-purpose models that treat MRI as one modality among many; large-scale, MRI-specific modelling and benchmarking remain limited, even though MRI’s low soft-tissue contrast leaves many boundaries effectively invisible on individual slices. We present SAMRI-3D, a benchmark and method for 3D MRI segmentation with SAM2. The SAMRI-3D benchmark is the largest MRI-only evaluation to date — 10,392 volumes from 34 datasets (27 public, 7 in-house) spanning 12 anatomical domains and 10+ sequences, with explicit seen/unseen splits. Freezing the image encoder and fine-tuning only the lightweight decoder and memory modules raises mean Dice from 0.58 (zero-shot SAM2) to 0.76, surpassing recent SAM-based medical models (SAMed-2 0.69, Medical-SAM2 0.49, SAM-Med3D 0.37) with strong statistical significance. To target invisible boundaries, we introduce Global Volume Tokens (GVT): persistent memory tokens trained with a Truncated Signed Distance Field (TSDF) reconstruction objective that is discarded at inference (zero added cost). This full model, SAMRI-3D, attains the best accuracy (0.78) and lowest variance across all 34 datasets and, uniquely, shows no drop on 8 held-out datasets (0.79 unseen vs. 0.78 seen); per-sequence analysis confirms the TSDF objective helps most where per-slice contrast is weakest. We will release the benchmark, code, and models in this paper.

Keywords: Medical image segmentation, MRI, SAM2, foundation model adaptation, 3D segmentation, benchmark.

1. Introduction

Magnetic resonance imaging (MRI) segmentation underpins diagnosis, treatment planning, and disease moni-

toring, yet it is difficult: intensities are non-uniform, contrast varies with protocol, and many important structures are small or weakly delineated. Manual delineation is accurate but costly and varies between experts. Supervised networks — whether 2D or 3D, from nnU-Net [13] and SwinUNETR [14] to specialised MRI models [15, 16, 17] — are accurate but must be retrained whenever the target or protocol changes. Promptable foundation models such as the Segment Anything Model (SAM) [1] offer an alternative — segmenting arbitrary targets from simple prompts without per-target retraining — but SAM is trained on natural images and transfers poorly to medical data.

Medical adaptations such as MedSAM [2] and the MRI-specific SAMRI [3] reduce this domain gap but remain 2D: applied to a volume they segment each slice independently, ignoring through-plane continuity. SAM2 [4] extends SAM to video, propagating predictions across frames through a streaming memory bank — a design that maps naturally onto 3D volumes via the “slice-as-frame” paradigm (Section 2.1). Still, SAM2 underperforms out-of-the-box on medical scans, and its finite memory bank loses information from distant slices in long volumes.

Existing SAM-for-3D adaptations [18, 19, 20, 21, 22, 23] instead add dedicated 3D modules or discard SAM2’s streaming memory entirely. We keep that memory and address two open problems. First, existing SAM2 medical evaluations each cover only a few organs or a single modality; no large-scale, MRI-specific benchmark exists. Second, SAM2’s memory bank provides only local inter-slice coherence and lacks global volumetric context — information from the whole volume that could guide ambiguous slices; recent memory improvements [7, 8] still operate within a local temporal window. We propose Global Volume Tokens (GVT) to inject whole-volume context into SAM2’s memory bank, trained with a Truncated Signed Distance Field (TSDF) objective that encodes boundary geometry rather than redundant binary masks.

Our contributions are: (1) the SAMRI-3D benchmark, the largest MRI-only SAM2 evaluation to date — 10,392

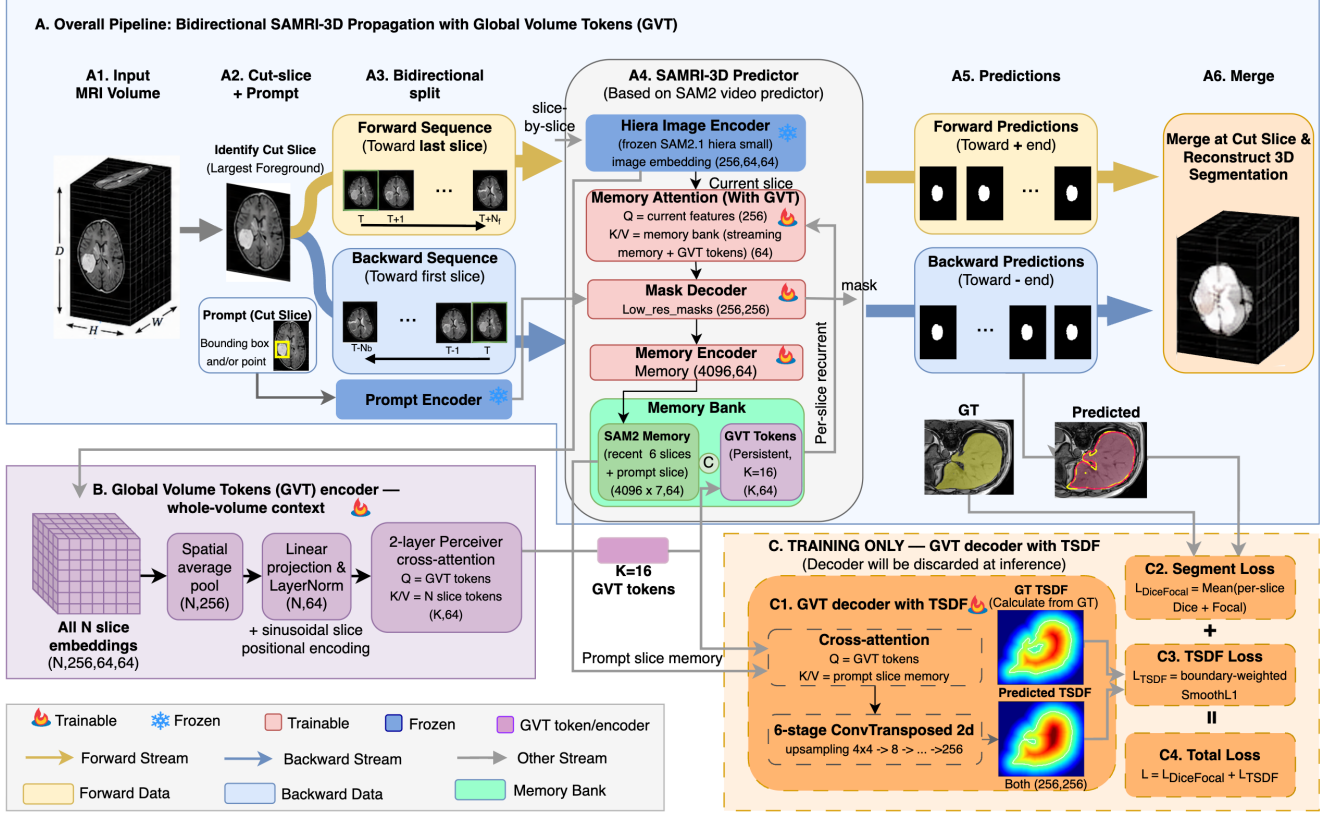


Figure 1. End-to-end SAMRI-3D workflow. (A) Given a 3D MRI volume and a box/point prompt on its cut slice (the slice with the largest target foreground), the model splits the volume into a forward and a backward sub-sequence and segments it slice by slice, propagating outward from the cut slice in both directions before merging the two halves into the final 3D mask. The predictor is SAM2 with a frozen image encoder and fine-tuned decoder and memory modules; at every step its memory attention reads a streaming memory bank that we augment with a persistent set of Global Volume Tokens (GVT). (B) The GVT encoder builds these tokens by compressing the embeddings of all slices into a compact summary of the whole volume, which is injected into the memory bank so that every slice—however far from the cut slice—sees global context. (C) During training only, a lightweight decoder uses the GVT tokens to reconstruct a Truncated Signed Distance Field (TSDF) of the target and is supervised against the ground-truth TSDF; this reconstruction objective forces the tokens to encode boundary geometry. The network is trained with the segmentation loss plus this TSDF loss, and the TSDF decoder is discarded at inference, adding no cost.

volumes from 34 datasets (27 public, 7 in-house), 12 anatomical domains, and 10+ sequences, with seen/unseen splits for zero-shot testing; (2) an efficient MRI adaptation of SAM2 that freezes the image encoder and fine-tunes only the lightweight decoder and memory modules, lifting mean Dice from 0.58 (SAM2 zero-shot) to 0.76 and already outperforming recent SAM-based medical models by statistically significant margins; (3) Global Volume Tokens (GVT) with a Truncated Signed Distance Field (TSDF) reconstruction objective — the first study of compact, persistent whole-volume memory tokens for SAM2; we compare training these tokens with segmentation loss alone, auxiliary binary reconstruction, and TSDF reconstruction, and show the TSDF variant is the most effective, giving the best accuracy (0.78) and lowest variance at zero inference overhead; and (4) robust generalisation and mechanistic analy-

sis — our model transfers to held-out datasets with no performance drop, and per-sequence analysis shows the reconstruction objective helps exactly where image boundaries are weakest.

2. Related Work

2.1. SAM2 and Its Use for 3D MRI

SAM2 [4] is a promptable segmentation model for both images and video, pairing a Hiera image encoder, a prompt encoder (points or boxes), and a mask decoder with a streaming memory bank. A video is segmented frame by frame: when predicting the current frame, a memory attention module conditions it on features stored from previously processed frames. Because adjacent frames are highly correlated, this lets a mask prompted on a few frames prop-

agate coherently through the whole sequence. The same property holds along the third axis of a 3D volume, where neighbouring slices are strongly correlated, so SAM2 fits 3D MRI under the slice-as-frame paradigm: the volume is treated as a video and a mask prompted on one slice is propagated slice by slice, with the memory bank carrying the inter-slice context that a purely 2D model would discard [5, 6]. This through-plane coherence is precisely why we build on SAM2 rather than a 2D segmenter; we detail the memory bank, and its limits, in Section 2.2. Medical adaptations of SAM2 fine-tune it or adjust its memory — MedSAM2 [9] on 455K+ image-mask pairs, SAMed-2 [8] with temporal adapters and confidence-driven memory, and RevSAM2 [10] with reverse propagation — but are evaluated predominantly on CT or limited MRI subsets.

2.2. Memory Bank and Global Context

The memory bank is the component we build on. A memory encoder stores each processed frame as a spatial feature map paired with its predicted mask, and the memory attention retrieves a fixed number of the most recent such entries when segmenting the current frame. Because SAM2 is built to segment unbounded video, this makes the bank a sliding window: it gives strong local coherence between neighbouring slices but evicts distant ones, so information from far parts of a long volume is lost. An MRI volume, however, is finite — at most a couple of hundred slices — so its entire content can be summarised once and made globally available, which is exactly what our Global Volume Tokens do (Section 3.3). Prior work instead improves memory quality within the local window — short/long-term banks (SLM-SAM2 [7]), self-sorting confidence-based selection (MedSAM-2 [11]), and confidence-driven replacement (SAMed-2 [8]) — but none provides whole-volume context. Compressing global context into a compact set of learnable tokens has precedent in MRI: Bran Lorenzana *et al.* [12] used global image tokens for transformer-based reconstruction.

2.3. Signed Distance Fields

A signed distance field (SDF) replaces a binary mask with each pixel’s signed distance to the nearest boundary — positive inside, negative outside, zero on the boundary — a dense, geometry-aware target that emphasises boundary localisation over per-pixel occupancy. Recent segmentation work exploits this through boundary-weighted losses (FocusSDF [24]), flow-matching (FlowSDF [25]), and SDF-supervised token compression (TokenSeg [26]). A full SDF is unbounded, so much of its loss is spent on voxels far from the boundary; a truncated SDF (TSDF) clamps the field to a narrow band around the boundary, concentrating supervision where it matters. We compute a per-slice TSDF from the ground-truth mask (see supplementary) and use it as the

reconstruction target for our Global Volume Tokens (Section 3.3).

3. Method

3.1. Problem Formulation

Given a 3D MRI volume $V \in \mathbb{R}^{D \times H \times W}$, we treat it as a D -frame sequence and identify the *cut slice* [6] — the slice with the largest target foreground — which gives the clearest view of the target and the most confident single-slice prediction, the best anchor for propagation. We split the volume at this slice into a forward and a backward sub-sequence (Fig. 1(A)), each processed independently by SAMRI-3D’s video predictor (fine-tuned from SAM2) with the cut slice as the prompted frame. During training the prompt is a bounding box derived from the ground-truth mask on the cut slice; at inference, boxes and/or points are supplied by the user or an upstream detector. Predictions propagate outward from the anchor in both directions and are merged at the cut slice into the full 3D mask, preserving the original NIfTI spatial metadata.

3.2. Architecture

Following the parameter-efficient strategy of our prior 2D adaptation SAMRI [3], we freeze the Hiera image encoder and the prompt encoder and fine-tune only the lightweight decoder and memory modules (mask decoder, memory encoder, and memory attention). Because the encoder is frozen, its embeddings are pre-computed once for the whole training set (~ 3.3 TB), avoiding repeated encoder passes and greatly reducing training cost. The model then operates on three complementary levels of spatial context: within-slice features from the frozen encoder; local inter-slice coherence from the streaming memory bank; and whole-volume context from the Global Volume Tokens (Section 3.3), which persist in the memory bank so that every slice, however far from the cut slice, sees a global summary.

3.3. Global Volume Tokens

To give every slice whole-volume context, we introduce Global Volume Tokens (GVT). The GVT encoder (Fig. 1(B)) takes the pre-computed Hiera features of all D slices, spatially average-pools each into a slice vector, adds a sinusoidal positional encoding along the volume axis, and compresses the sequence with a 2-layer Perceiver-style [29] cross-attention into $K = 16$ persistent tokens (dimension $C_{\text{mem}} = 64$). These tokens are appended to the memory bank at every propagation step, so — unlike the sliding-window memory — they remain visible however far a slice is from the cut slice.

Trained through the segmentation loss from SAM2 alone, however, the tokens tend to become redundant with

the memory bank: the memory attention’s softmax competition lets the per-frame entries dominate, leaving the tokens little direct signal. We therefore supervise them with an auxiliary reconstruction objective (Fig. 1(C)). A training-only decoder cross-attends the GVT tokens (Q) to the prompted-slice memory features (K/V) — a natural label identifier, so the same architecture generalises across all structures without a label-embedding table — and upsamples them, through 6-stage ConvTranspose2d from 4×4 to 256×256 , into a Truncated Signed Distance Field (TSDF) of the target. This is compared to the ground-truth TSDF (truncation radius $\tau = 32$ px) by a boundary-weighted SmoothL1 loss that upweights voxels near the boundary shell. Because a TSDF encodes boundary geometry at every voxel rather than binary occupancy, it forces the tokens to carry shape information complementary to the appearance-based memory — unlike a binary reconstruction target, which duplicates what the mask decoder already produces. The decoder is discarded at inference, so GVT adds no test-time cost; the supplementary gives full layer configurations, tensor shapes, and the trainable/frozen breakdown. The model is trained with

$$\mathcal{L} = \mathcal{L}_{\text{DiceFocal}} + \mathcal{L}_{\text{TSDF}}. \quad (1)$$

4. Experiments and Results

4.1. Dataset

The SAMRI-3D benchmark comprises 10,392 3D MRI volumes from 34 datasets (27 public, 7 in-house) spanning 12 anatomical domains. By comparison, prior SAM2 medical studies are evaluated on far fewer datasets — for example, a single knee-MRI dataset [5], or BraTS plus four Medical Segmentation Decathlon tasks [6] — so no existing study assesses SAM2 on MRI at this scale. We partition the 34 datasets into 26 seen and 8 unseen datasets; each dataset is split 80/10/10 into train/validation/test. Models are trained and validated only on the seen datasets’ train/validation subsets and evaluated on the test subsets of both the seen and unseen datasets, the latter measuring zero-shot generalisation to datasets unseen during training. Table 1 summarises the benchmark by anatomical domain; the complete list of all 34 datasets with per-dataset sequence, volume count, and seen/unseen split is provided in Table S5 of the supplementary material.

4.2. Implementation Details

All models use the SAM2.1 Hiera-Small backbone with pre-computed slice embeddings, trained from these embeddings with AdamW [28] (learning rate 5×10^{-6} , weight decay 0.01). Each volume is trained with an 80%/20% mix of volumetric sequences and 2D slice batches — the 2D batches preserve strong per-slice segmentation, which sup-

Domain	# DS	# Vol.	Sequences
Brain	8	5,520	FLAIR, T1, T1CE, ADC, DWI, T2
Prostate	7	2,234	ADC, T2
Knee/MSK	8	1,053	FLASH, PD, T2, DESS, TSE
Cardiac	3	340	Cine, LGE
Whole-body	1	263	Mixed
Hippocampus	1	260	T1
Vestibular	1	227	T1CE
Hip	1	211	T2
Spine	1	172	T2
Abdominal	1	60	T1
Shoulder	1	45	T2
Gynecology	1	7	T2
Total	34	10,392	10+ MRI sequences

Table 1. SAMRI-3D Benchmark Summary by Anatomical Domain.

plies the reliable prompted-slice predictions that condition subsequent propagation.

To study how the GVT tokens should be trained, we compare four models (Fig. S1 in the supplementary): a 2D-slice baseline without GVT (GVT-baseline); GVT tokens trained through the segmentation loss alone (GVT-v1); with an auxiliary low-resolution binary reconstruction (GVT-v2); and with the TSDF reconstruction (GVT-v3, our SAMRI-3D). All use the per-slice Dice + Focal loss [27]; GVT-v3 adds the TSDF term of Eq. (1).

We report mean 3D Dice ($\pm$ std across the 34 datasets), Hausdorff Distance (HD), and Mean Surface Distance (MSD); per dataset we take the median across volumes, then the mean and std of these per-dataset medians. Significance uses two-sided paired Wilcoxon signed-rank tests at the per-sample ($n \approx 3,376$), per-dataset ($n = 34$), and per-label ($n = 49$) levels; full statistics in Table S12.

4.3. Main Results

Table 2 reports the headline comparison: three prior SAM-based medical methods, zero-shot SAM2-Small, and our SAMRI-3D. Because SAMRI-3D is built on the Hiera-Small backbone, SAM2-Small is its backbone-matched zero-shot reference. SAMRI-3D attains the best mean accuracy and lowest variance on every metric, lifting mean Dice from 0.58 (zero-shot SAM2-Small) to 0.78 and clearly surpassing SAMed-2, Medical-SAM2, and SAM-Med3D. All improvements over the comparison methods and SAM2-Small are significant at the per-sample and per-label levels (full statistics in Table S12). Two further SAM2-based adaptations produced near-zero Dice on this benchmark and are omitted. Backbone size has little effect (Section 4.4) and the token-training objective is analysed in Section 4.5. Per-dataset and per-target breakdowns are in Tables S6–S11.

4.4. SAM2 Scale Analysis

The aim of this analysis is to justify our backbone choice. SAM2’s encoder scale has little effect on zero-shot MRI accuracy: across the four SAM2 sizes (Table 3), a $5.8 \times$

Method	Dice $\uparrow$	HD $\downarrow$	MSD $\downarrow$
SAM-Med3D [21]	0.37 ± 0.21	29.04 ± 16.70	6.39 ± 2.78
Medical-SAM2 [11]	0.49 ± 0.22	33.51 ± 21.81	5.92 ± 4.26
SAMed-2 [8]	0.69 ± 0.15	22.28 ± 18.82	2.91 ± 1.62
SAM2-Small	0.58 ± 0.24	28.24 ± 20.42	6.32 ± 4.46
SAMRI-3D (Ours)	0.78 ± 0.14	13.85 ± 9.01	2.11 ± 0.96

Table 2. Main comparison: three prior SAM-based medical methods, the backbone-matched zero-shot SAM2-Small, and SAMRI-3D (mean $\pm$ std across 34 datasets). SAM2 backbone sizes and GVT ablations are in Tables 3 and 4.

Backbone	Params	Dice $\uparrow$	HD $\downarrow$	MSD $\downarrow$
SAM2-Tiny	38.9M	0.57 ± 0.23	27.79 ± 20.19	6.73 ± 5.07
SAM2-Small	46.0M	0.58 ± 0.24	28.24 ± 20.42	6.32 ± 4.46
SAM2-Base+	80.8M	0.56 ± 0.24	29.95 ± 21.39	6.23 ± 4.22
SAM2-Large	224.4M	0.56 ± 0.24	30.22 ± 20.87	7.83 ± 5.43

Table 3. Zero-shot SAM2 across backbone sizes. Accuracy is flat despite a $5.8\times$ parameter increase, so we fine-tune from SAM2-Small.

increase in parameters (38.9M to 224.4M) brings no improvement — all four variants fall within a narrow 0.56–0.58 mean Dice band, and the two largest models are in fact marginally worse, indicating that the extra capacity trained on natural video transfers poorly to MRI. We therefore avoid a large backbone and fine-tune from SAM2-Small, whose zero-shot accuracy is the best (0.58) at near-minimal parameter cost (46.0M, versus 224.4M for Large), giving the best accuracy–efficiency trade-off for fine-tuning. This size-insensitivity also supports our frozen-encoder design: since a larger encoder does not help, we freeze the Hiera encoder and adapt only the decoder and memory modules.

4.5. GVT Ablation Analysis

Table 4 compares the four training strategies. Adding GVT tokens without a reconstruction objective (GVT-v1) or with a binary reconstruction (GVT-v2) gives little or only marginal gain over the baseline, whereas the TSDF reconstruction (GVT-v3, our SAMRI-3D) improves every metric (Dice 0.78 ± 0.14 , HD 13.85 ± 9.01 , MSD 2.11 ± 0.96) and yields the lowest Dice variance (± 0.14 vs. ± 0.17 for the baseline). The gain over the strong baseline is small but significant on Dice at the per-sample ($p = 9.8 \times 10^{-14}$) and per-label ($p = 4.1 \times 10^{-3}$, 71% win-rate) levels, though not at the coarser per-dataset level ($p = 0.22$); GVT-v3 also exceeds GVT-v1 at all levels ($p \leq 1.1 \times 10^{-2}$) while being statistically indistinguishable from GVT-v2, its edge lying in mean accuracy and reduced variance. This progression is exactly what the design in Section 3.3 predicts: tokens trained by the segmentation loss alone stay redundant with the memory bank, a binary reconstruction target duplicates what the mask decoder already learns, and only the TSDF’s boundary geometry supplies information genuinely complementary to SAM2’s appearance-based memory.

Method	Objective	Dice $\uparrow$	HD $\downarrow$	MSD $\downarrow$
GVT-baseline	—	0.76 ± 0.17	14.21 ± 12.32	2.39 ± 1.70
GVT-v1	seg	0.76 ± 0.16	14.46 ± 10.09	2.33 ± 1.15
GVT-v2	binary	0.77 ± 0.16	13.94 ± 10.34	2.25 ± 1.40
GVT-v3 (SAMRI-3D)	TSDF	0.78 ± 0.14	13.85 ± 9.01	2.11 ± 0.96

Table 4. GVT ablation: how the global tokens are trained, each objective additive on the per-slice segmentation loss. Objective — seg: segmentation loss only; binary: + binary-mask reconstruction; TSDF: + TSDF reconstruction.

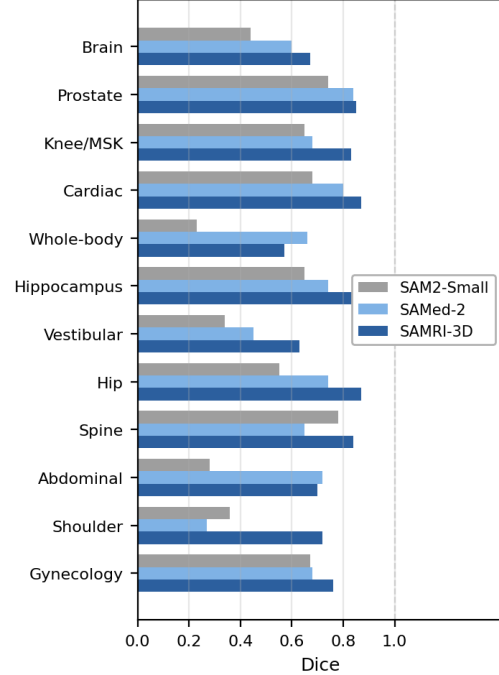


Figure 2. Per-domain Dice (mean of per-dataset medians) for SAMRI-3D versus zero-shot SAM2-Small and the strongest external method (SAMed-2) across the 12 benchmark domains. SAMRI-3D leads in 10 of 12, trailing SAMed-2 only on the whole-body and abdominal multi-organ cases.

4.6. Per-Domain Highlights

SAMRI-3D leads in 10 of 12 anatomical domains, trailing SAMed-2 only on the whole-body and abdominal multi-organ cases (Fig. 2, Table S3). It is strongest on well-defined structures and, more importantly, rescues low-contrast targets that zero-shot SAM2 largely misses — white-matter hyperintensities improve from 0.09 to 0.54 and the hippocampus from 0.65 to 0.84 — while a few small, diffuse structures (scapula cartilage, non-enhancing tumour) remain hard for every method.

4.7. Per-MRI-Sequence Analysis

Because its TSDF objective supplies boundary geometry, SAMRI-3D should help most where image evidence is weakest, so its effect should track imaging contrast. Fig. 3 confirms this: across 11 single-sequence categories, gains

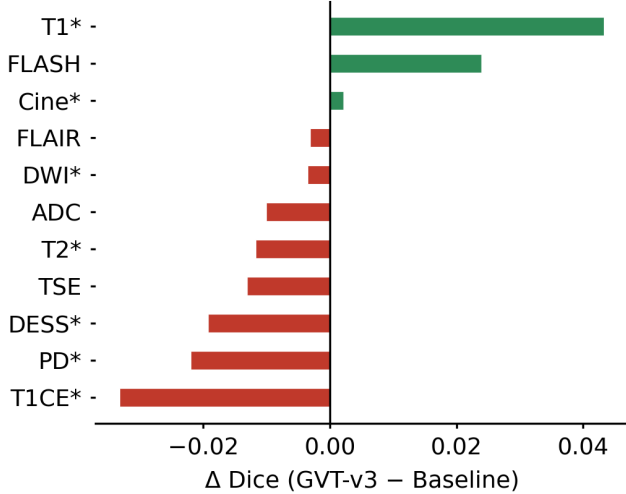


Figure 3. Per-MRI-sequence improvement of SAMRI-3D over the baseline (median Dice difference). Green: improvement; red: regression; * marks $p < 0.05$ (paired Wilcoxon). Single-sequence subset (2,156 samples).

concentrate on low-contrast sequences with diffuse boundaries (T1, FLASH; pooled $+0.035$ mean Dice, $p = 1.4 \times 10^{-6}$) and are neutral-to-slightly-negative on high-contrast sequences whose boundaries are already sharp (DESS, T1CE). This sequence dependence is direct evidence for the invisible-boundary motivation behind the TSDF objective. Full per-sequence values are reported in Table S4.

4.8. Zero-Shot Generalisation to Unseen Datasets

Because the 8 unseen datasets contribute no training data, their test subsets probe generalisation to datasets absent from training. Fig. 4 compares seen- and unseen-dataset performance for all models. SAMRI-3D attains the highest unseen Dice (0.79) and is the only strong model with no generalisation gap ($+0.01$ relative to seen data), whereas SAMed-2 and the fine-tuned baseline both degrade on unseen data. On the unseen datasets it significantly exceeds SAMed-2 and zero-shot SAM2 ($p < 0.02$); its edge over the fine-tuned baseline is positive but not significant at $n = 8$. Gains on unseen musculoskeletal data (e.g., MSK-FLASH, $0.38 \rightarrow 0.71$) show the TSDF objective transfers robustly.

4.9. Qualitative Results

Fig. 5 compares SAMRI-3D against zero-shot SAM2 (Small), SAM-Med3D, Medical-SAM2, and SAMed-2. Successful cases are shown in triplanar (axial, coronal, sagittal) views for liver, spleen, and cartilage; SAMRI-3D produces the most complete and boundary-accurate masks, recovering structure that the zero-shot and external models under-segment or miss. The failure cases, shown as single slices, illustrate that all methods, including SAMRI-

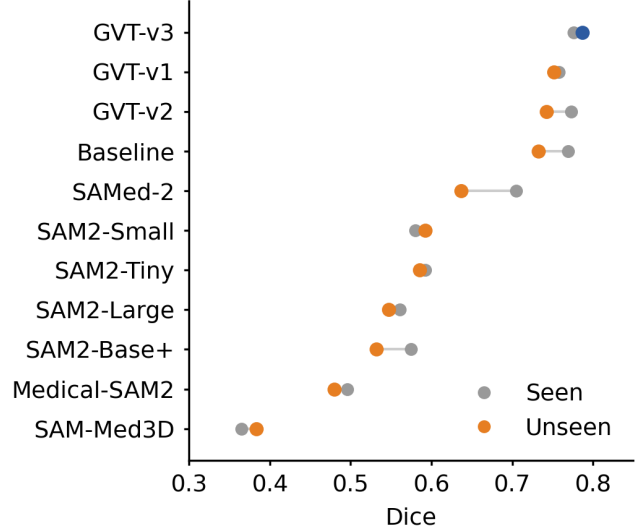


Figure 4. Zero-shot generalisation: per-model seen vs. unseen Dice (mean across datasets). SAMRI-3D (highlighted) shows no generalisation gap. Seen: 26 datasets used for training/validation; unseen: 8 datasets excluded from training.

3D, struggle on extremely small or diffuse targets, such as a scattered brain tumour and a thin cartilage, where every model attains low and closely clustered Dice, reflecting genuine task difficulty rather than a model-specific weakness. Additional qualitative examples across more anatomical domains are provided in Fig. S2 of the supplementary material.

5. Discussion

Frozen-encoder fine-tuning works because of the large domain gap between natural video and MRI: the frozen Hierarchical encoder supplies general visual features while the trained decoder and memory adapt them to medical imaging, preserving the encoder’s generalisation.

The finding that model size does not matter for zero-shot medical performance is notable. All SAM2 variants, from Tiny to Large, achieve similarly low zero-shot accuracy, suggesting that the bottleneck is not how well the model “sees” the image (the frozen encoder already extracts strong general features) but how well the decoder and memory interpret those features for the medical domain. Among existing SAM-based medical methods, SAMed-2’s relatively strong performance confirms the value of adapter-based fine-tuning.

The GVT progression also clarifies how global context should enter SAM2’s memory: only the TSDF reconstruction, which encodes boundary geometry rather than appearance, adds information that survives the memory attention’s softmax competition (Section 4.5). Its improvement is modest but consistent across all metrics and datasets, with notably reduced variance; as the significance tests show, the

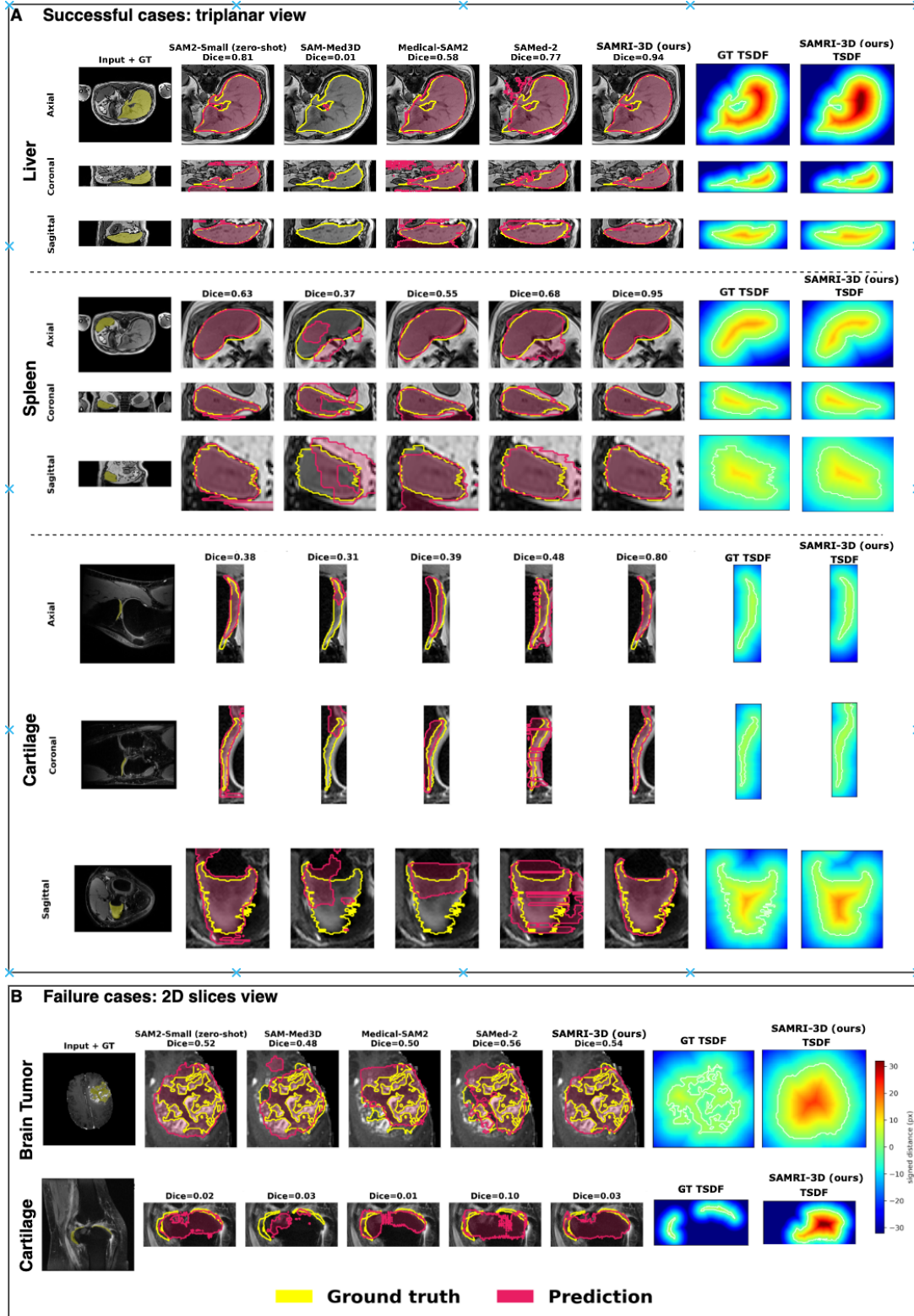


Figure 5. Qualitative segmentation results. Top rows A: representative successful cases shown in triplanar (axial, coronal, sagittal) views, with ground-truth and predicted masks overlaid; the two rightmost columns show the target’s ground-truth TSDF and the TSDF reconstructed by SAMRI-3D’s training-only decoder. Bottom rows B: two characteristic failure modes — (i) a brain tumour whose ground truth is highly fragmented, which almost all models instead predict as a single contiguous region; and (ii) a thin cartilage whose bounding-box prompt spans a large area, producing an over-segmentation on the prompted slice that then propagates into an incorrect 3D volume.

gain is robust per-sample and per-label but not per-dataset, so we position SAMRI-3D as the lowest-variance, best-mean model rather than claiming a large dataset-level improvement.

Several directions remain open. Gated attention that bypasses softmax competition could amplify SAMRI-3D's gains by modulating features directly rather than competing for attention weight, and adapting the frozen encoder may lift its performance ceiling. Small-lesion segmentation is still hard (e.g., adrenal Dice ≤ 0.10), and comparison with fully-supervised methods such as nnU-Net would further contextualise the results. Finally, the two failure modes in Fig. 5 — fragmented ground truth collapsed into one contiguous mask, and over-segmentation from over-sized box prompts propagating through the memory bank — motivate scatter-aware prompting and tighter prompt-region handling.

6. Conclusion

We present SAMRI-3D, a comprehensive benchmark and method for adapting SAM2 to 3D MRI segmentation. Our best model, SAMRI-3D with TSDF reconstruction, achieves the best mean accuracy with the lowest variance across 34 diverse datasets, substantially outperforming three existing SAM-based medical segmentation methods. We show that TSDF-guided Global Volume Tokens provide consistent improvements where simpler token formulations fall short, pointing toward structured geometric priors as a principled approach for global context integration. The SAMRI-3D benchmark, comprising 10,392 volumes across 12 anatomical domains, provides a standardised evaluation framework for future MRI segmentation research.

Code Availability

Code and pre-trained models are publicly available at <https://github.com/wangzhaomxy/SAMRI-3D>.

Acknowledgements

This work was partially supported by a National Health and Medical Research Council (NHMRC) Ideas Grant (Award No. APP2001734). During manuscript preparation the authors used a large language model to check grammar and refine language; all such content was reviewed and edited by the authors, who take full responsibility for it.

References

- [1] A. Kirillov *et al.*, “Segment anything,” in *Proc. ICCV*, 2023. **1**
- [2] J. Ma *et al.*, “Segment anything in medical images,” *Nature Communications*, vol. 15, 2024. **1**
- [3] Z. Wang *et al.*, “SAMRI: Segment any MRI,” arXiv:2510.26635, 2025. **1, 3**
- [4] N. Ravi *et al.*, “SAM 2: Segment anything in images and videos,” arXiv:2408.00714, 2024. **1, 2**
- [5] A. S. Yu *et al.*, “Novel adaptation of video segmentation to 3D MRI: Efficient zero-shot knee segmentation with SAM2,” arXiv:2408.04762, 2024. **3, 4**
- [6] C. Shen *et al.*, “Interactive 3D medical image segmentation with SAM 2,” arXiv:2408.02635, 2024. **3, 4**
- [7] Y. Chen *et al.*, “SLM-SAM 2: Short-long memory SAM 2 for medical image segmentation,” arXiv:2505.01854, 2025. **1, 3**
- [8] Z. Cheng *et al.*, “SAMed-2: Adapting SAM2 for medical image segmentation,” in *Proc. MICCAI*, 2025. **1, 3, 5**
- [9] J. Ma *et al.*, “Segment anything in medical images and videos: Benchmark with application to echocardiography segmentation,” arXiv, 2025. **3**
- [10] Y. Bai *et al.*, “RevSAM2: Reverse propagation for medical image segmentation,” arXiv:2409.04298, 2024. **3**
- [11] J. Zhu *et al.*, “Medical SAM 2: Segment medical images as video via segment anything model 2,” arXiv:2408.00874, 2024. **3, 5**
- [12] M. Bran Lorenzana *et al.*, “Transformer compressed sensing via global image tokens,” in *Proc. IEEE ICIP*, 2022, pp. 3011–3015. **3**
- [13] F. Isensee *et al.*, “nnU-Net: a self-configuring method for deep learning-based biomedical image segmentation,” *Nature Methods*, vol. 18, pp. 203–211, 2021. **1**
- [14] A. Hatamizadeh *et al.*, “Swin UNETR: Swin transformers for semantic segmentation of brain tumours in MRI images,” in *BrainLes, MICCAI*, 2022. **1**
- [15] W. Dai *et al.*, “CAN3D: Fast 3D medical image segmentation via compact context aggregation,” *Medical Image Analysis*, vol. 82, art. 102562, 2022. **1**
- [16] B. Woo *et al.*, “Semantic segmentation of 3D medical images through a kaleidoscope: Data from the Osteoarthritis Initiative,” in *Proc. MIDL, PMLR*, vol. 227, 2023, pp. 905–917. **1**
- [17] B. Woo *et al.*, “Automated anomaly-aware 3D segmentation of bones and cartilages in knee MR images from the Osteoarthritis Initiative,” *Medical Image Analysis*, vol. 93, art. 103089, 2024. **1**
- [18] C. Chen *et al.*, “MA-SAM: Modality-agnostic SAM adaptation for 3D medical image segmentation,” *Medical Image Analysis*, vol. 98, 2024. **1**
- [19] S. Gong *et al.*, “3DSAM-adapter: Holistic adaptation of SAM from 2D to 3D for promptable medical image segmentation,” *Medical Image Analysis*, vol. 98, 2024. **1**
- [20] J. Wu *et al.*, “Medical SAM adapter: Adapting segment anything model for medical image segmentation,” *Medical Image Analysis*, vol. 102, 2024. **1**
- [21] S. Wang *et al.*, “SAM-Med3D: Towards general-purpose segmentation models for volumetric medical images,” arXiv, 2024. **1, 5**
- [22] A. H. Abdusalamov *et al.*, “3DMedSAM: Cross-dimensional adaptation of SAM for volumetric medical image segmentation,” *Biocybernetics and Biomedical Eng.*, 2024. **1**
- [23] H. Xiong *et al.*, “SAM2-UNet: Segment anything 2 makes strong encoder for natural and medical image segmentation,” arXiv:2408.08870, 2024. **1**
- [24] U. Shafique *et al.*, “FocusSDF: Boundary-aware learning via signed distance supervision,” arXiv:2511.11864, 2025. **3**
- [25] F. Bogensperger *et al.*, “FlowSDF: Flow matching for medical image segmentation using distance transforms,” *Int. J. Computer Vision*, 2025. **3**

- [26] Z. Zeng *et al.*, “TokenSeg: Efficient 3D medical image segmentation via hierarchical visual token compression,” *Int. J. Computer Vision*, 2025. 3
- [27] T.-Y. Lin *et al.*, “Focal loss for dense object detection,” in *Proc. ICCV*, 2017. 4
- [28] I. Loshchilov and F. Hutter, “Decoupled weight decay regularization,” in *Proc. ICLR*, 2019. 4
- [29] A. Jaegle *et al.*, “Perceiver: General perception with iterative attention,” in *Proc. ICML*, 2021. 3

Supplementary Material

A. TSDF Computation

For each slice we convert the ground-truth binary mask into a Truncated Signed Distance Field (TSDF). Let d_{in} and d_{out} be the Euclidean distance transforms of the foreground and background. The signed distance field

$$\text{SDF} = d_{\text{in}} - d_{\text{out}} \quad (2)$$

is positive inside the structure, negative outside, and zero on the boundary. We truncate it to a narrow band,

$$\text{TSDF} = \text{clip}(\text{SDF}, -\tau, +\tau), \quad \tau = 32 \text{ px}, \quad (3)$$

so that an all-background slice maps to the constant $-\tau$; targets are computed at 256×256 resolution. The TSDF decoder is supervised with a boundary-weighted Smooth-L1 loss

$$\mathcal{L}_{\text{TSDF}} = \text{mean}(w \odot \text{SmoothL1}(\hat{y}, y)), \quad w = 1 + \alpha \exp\left(-\frac{y^2}{2\sigma^2}\right), \quad (4)$$

where the weight w up-weights voxels near the boundary shell ($|y| \approx 0$), with $\alpha = 4$ and $\sigma = \tau/3 \approx 10.7$. For multi-label volumes the per-label loss is divided by $\sqrt{\text{area}}$ to balance small and large structures.


B. Architecture and Implementation Details

This section specifies the SAMRI-3D components and their tensor shapes for reproducibility; all values are for the SAM2.1 Hiera-Small backbone. Below, N is the number of slices in a volume, K the number of Global Volume Tokens, C_{mem} the SAM2 memory dimension, and C_{hid} the encoder/hidden dimension.

Key dimensions. The constants used throughout are summarised below.


Symbol	Meaning	Value
N	slices per volume	50–200
K	Global Volume Tokens	16
C_{mem}	SAM2 memory dimension	64
C_{hid}	encoder / hidden dimension	256
—	input resolution	1024×1024
—	encoder feature map	64×64
M	memory tokens per frame (64×64)	4096
—	TSDF resolution	256×256
τ	TSDF truncation radius	32 px
num_maskmem	memory sliding window	≈ 7 frames
max_slices	per-video slice cap (training)	32

GVT encoder. The encoder compresses all N pre-computed Hiera slice embeddings into K global tokens, run once per volume. Each slice embedding ($256 \times 64 \times 64$) is spatially average-pooled to a 256-d vector; a linear projection maps it to 64-d and a fixed sinusoidal positional encoding along the slice axis is added, giving an $(N, 64)$ sequence. $K = 16$ learnable query tokens then cross-attend to this sequence through a 2-layer Perceiver-style Transformer decoder (4 heads, feed-forward width 256), producing K global tokens of dimension 64. The encoder has $\approx 120\text{K}$ parameters.

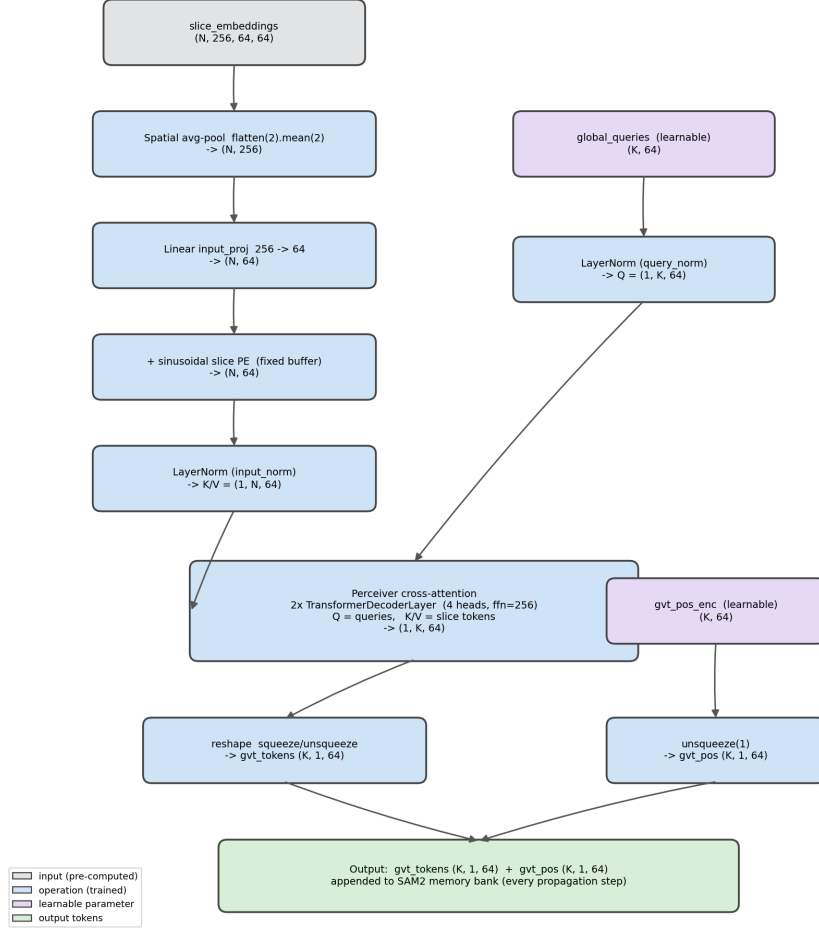


GVT encoder data flow with tensor shapes.

Memory-bank fusion. SAM2s memory attention is wrapped so that, at every propagation step, the cached GVT tokens are concatenated onto the memory bank: the memory sequence grows from $(M_{\text{total}}, B, 64)$ to $(M_{\text{total}}+K, B, 64)$. Here M_{total} is SAM2s sliding window of ≈ 7 recent frames (each $64 \times 64 = 4096$ tokens) plus object-pointer tokens; the GVT tokens are treated like object-pointer tokens and excluded from spatial rotary position encoding. The current-frame query and attention output are $C_{\text{hid}}=256$ -dimensional, whereas the memory keys/values (and the GVT tokens) are $C_{\text{mem}}=64$ -dimensional; cross-attention projects $64 \rightarrow 256$ internally. Because the same $K=16$ tokens are appended for every slice, even the most distant slice attends to an identical whole-volume summary.

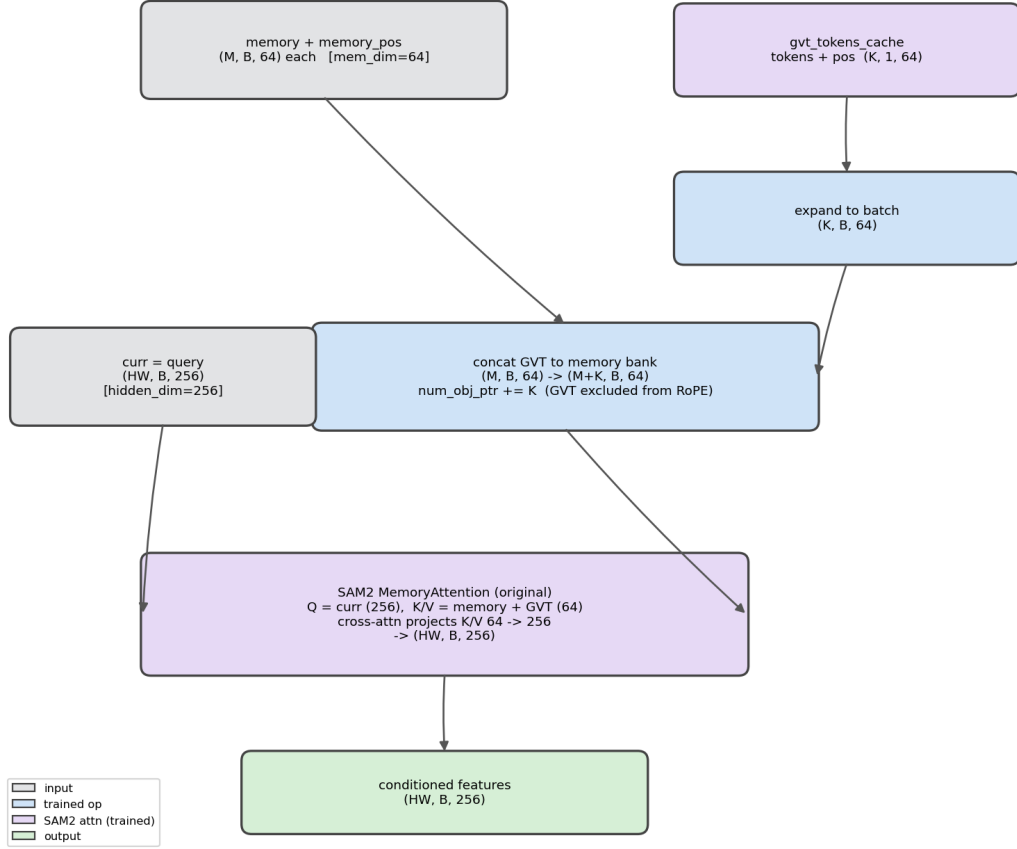

Memory-bank fusion: the K tokens are concatenated onto SAM2s memory bank before the original memory attention.

TSDF decoder (training only). For each label and slice, the K GVT tokens (each given an additive learned slice positional embedding) cross-attend to the prompted-slice memory features (4096×64). The conditioned tokens are flattened ($16 \times 64 = 1024$) and linearly projected to a $64 \times 4 \times 4$ spatial seed, then upsampled by six ConvTranspose2d stages (channels $64 \rightarrow 64 \rightarrow 64 \rightarrow 32 \rightarrow 32 \rightarrow 16 \rightarrow 8$; spatial $4 \rightarrow 256$, each with GroupNorm and GELU) and a final 1×1 convolution to a single-channel 256×256 TSDF (linear output). It has $\approx 30K$ parameters and is discarded at inference.

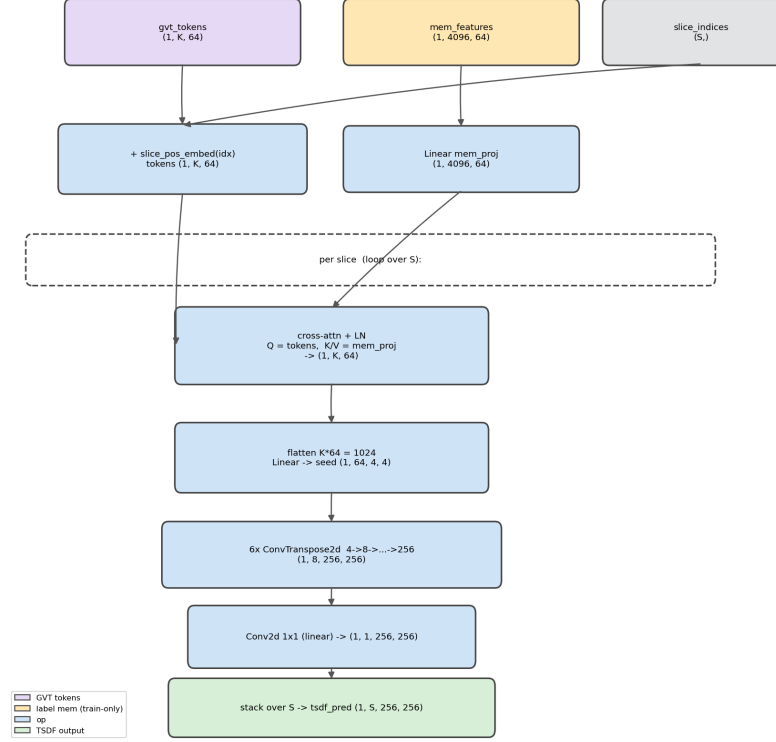

TSDF decoder (training only): per-slice cross-attention conditioning then a six-stage $4 \rightarrow 8 \rightarrow \dots \rightarrow 256$ ConvTranspose stack.

Trainable and frozen components. The image and prompt encoders are frozen; only the lightweight decoder, memory, and GVT modules are trained.

Component	Status	Params
Hiera image encoder	Frozen	—
Prompt encoder	Frozen	—
Mask decoder	Trained	SAM2
Memory encoder	Trained	SAM2
Memory attention	Trained	SAM2
GVT encoder	Trained	$\approx 120\text{K}$
TSDF decoder	Trained; discarded at inference	$\approx 30\text{K}$

Training. Models are trained from the pre-computed float16 embeddings with AdamW (learning rate 5×10^{-6} , weight decay 0.01), batch size 2, and AMP mixed precision, for 50 epochs on three GPUs. Each volume uses an 80%/20% mix of volumetric sequences and 2D slice batches, with video processing capped at $\text{max_slices} = 32$ per sub-sequence. The TSDF objective uses $\tau = 32$ px, $\alpha = 4$, and $\sigma = \tau/3 \approx 10.7$.

C. GVT Ablation Strategies

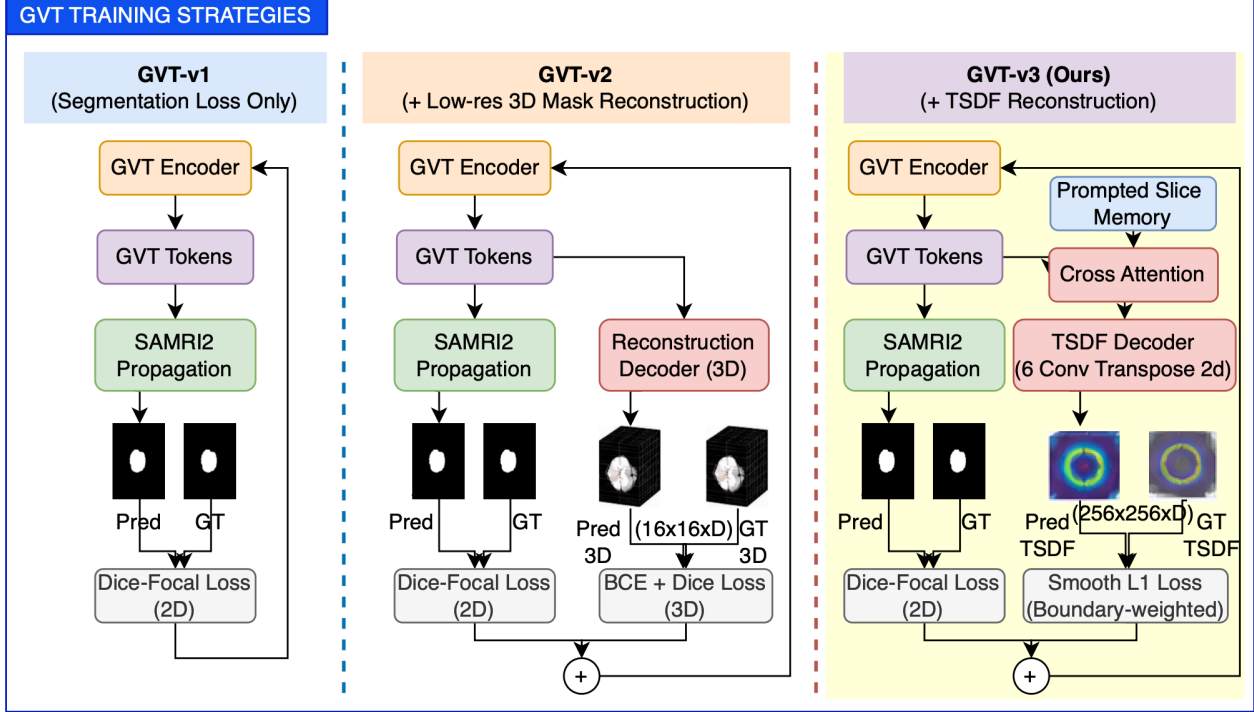


Figure S1. GVT ablation strategies. GVT-v1 trains the tokens through the SAM2 segmentation loss only; GVT-v2 adds a low-resolution binary-mask reconstruction; GVT-v3 replaces it with a Truncated Signed Distance Field (TSDF) reconstruction. The auxiliary decoder is discarded at inference.

D. Model Scaling and Efficiency

Table S1. SAM2 zero-shot performance by model size.

Model	Params	Mean Dice
sam2_tiny	38.9M	0.571
sam2_small	46.0M	0.583
sam2_b+	80.8M	0.564
sam2_large	224.4M	0.558

Table S2. Efficiency comparison of GVT variants (parameter overhead).

Variant	GVT Params	Aux Decoder	Total Overhead
Baseline	—	—	0
GVT-v1	120K	—	<0.9%
GVT-v2	120K	85K	<1.5%
SAMR1-3D	120K	30K	<1.1%

E. Per-Domain and Per-Sequence

Table S3. Per-domain Dice (mean $\pm$ SD of per-dataset median Dice across the 12 anatomical domains).

Domain	n	SAM2-Small	SAMed-2	SAMRI-3D
Brain	8	0.44 ± 0.29	0.60 ± 0.16	0.67 ± 0.17
Prostate	7	0.74 ± 0.15	0.84 ± 0.10	0.85 ± 0.11
Knee/MSK	8	0.65 ± 0.23	0.68 ± 0.05	0.83 ± 0.09
Cardiac	3	0.68 ± 0.21	0.80 ± 0.02	0.87 ± 0.02
Whole-body	1	0.23 ± 0.15	0.66 ± 0.16	0.57 ± 0.13
Hippocampus	1	0.65 ± 0.13	0.74 ± 0.04	0.84 ± 0.05
Vestibular	1	0.34 ± 0.11	0.45 ± 0.05	0.63 ± 0.09
Hip	1	0.55 ± 0.21	0.74 ± 0.03	0.87 ± 0.04
Spine	1	0.78 ± 0.07	0.65 ± 0.02	0.84 ± 0.02
Abdominal	1	0.28 ± 0.10	0.72 ± 0.04	0.70 ± 0.10
Shoulder	1	0.36 ± 0.16	0.27 ± 0.06	0.72 ± 0.18
Gynecology	1	0.67 ± 0.00	0.68 ± 0.00	0.76 ± 0.00

Table S4. Per-sequence Dice (mean ± SD of per-dataset median Dice across datasets with each sequence).

Sequence	n	SAM2-Small	SAMed-2	SAMRI-3D-2dloss	SAMRI-3D
T1	4	0.43 ± 0.30	0.65 ± 0.18	0.71 ± 0.16	0.72 ± 0.13
T1CE	2	0.51 ± 0.23	0.46 ± 0.01	0.71 ± 0.16	0.72 ± 0.12
FLAIR	3	0.39 ± 0.37	0.62 ± 0.20	0.60 ± 0.28	0.61 ± 0.25
DWI	1	0.58 ± 0.22	0.73 ± 0.06	0.81 ± 0.15	0.80 ± 0.13
ADC	3	0.56 ± 0.26	0.77 ± 0.15	0.72 ± 0.21	0.74 ± 0.17
T2	10	0.70 ± 0.17	0.72 ± 0.19	0.82 ± 0.12	0.84 ± 0.09
Cine	1	0.85 ± 0.07	0.79 ± 0.06	0.85 ± 0.07	0.85 ± 0.06
DESS	3	0.75 ± 0.04	0.67 ± 0.05	0.90 ± 0.02	0.89 ± 0.03
TSE	1	0.55 ± 0.08	0.59 ± 0.03	0.72 ± 0.08	0.70 ± 0.08
FLASH	2	0.38 ± 0.37	0.71 ± 0.04	0.58 ± 0.28	0.75 ± 0.06
PD	1	0.79 ± 0.11	0.71 ± 0.00	0.90 ± 0.01	0.88 ± 0.02

F. Benchmark Datasets

Table S5. Complete list of the 34 SAMRI-3D benchmark datasets.

Domain	Dataset	Sequence	# Vol.	# Labels	Split
Brain	MSD-BrainTumour [1]	FLAIR	1,936	3	Seen
Brain	BraTS-FLAIR [2]	FLAIR	954	1	Seen
Brain	BraTS-T1 [2]	T1w	954	1	Seen
Brain	BraTS-T1CE [2]	T1CE	890	1	Seen
Brain	ISLES2022-DWI [3]	DWI	224	1	Seen
Brain	ISLES2022-ADC [3]	ADC	224	1	Seen
Brain	WMH-T1 [4]	T1w	169	3	Seen
Brain	WMH-FLAIR [4]	FLAIR	169	3	Seen
Prostate	PI-CAI [5]	T2w	1,315	1	Seen
Prostate	Prostate-T2 [6]	T2w	338	1	Seen
Prostate	Prostate-ADC [7]	ADC	285	1	Seen
Prostate	QIN-Prostate [8]	T2w	90	1	Seen
Prostate	PROMISE12 [9]	T2w	80	1	Seen
Prostate	MSD-Prostate [1]	ADC	64	2	Seen
Prostate	QIN-Lesion [8]	T2w	62	1	Seen
Knee	ZIB-OAI [10]	DESS	507	4	Seen
Knee	OAI-iMorphics-DESS †	DESS	174	6	Unseen
Knee	OAI-iMorphics-TSE †	TSE	174	6	Unseen
Knee	OAI-iMorphics-FLASH †	FLASH	88	6	Unseen
Knee	OAI-AKOA [11]	DESS	38	6	Seen
Knee	MSK-T2 †	T2w	31	6	Seen
Knee	MSK-FLASH †	FLASH	30	6	Unseen
Knee	MSK-PD †	PD	11	6	Unseen
Cardiac	ACDC [12]	Cine bSSFP	300	3	Unseen
Cardiac	MSD-Heart [1]	LGE	20	1	Seen
Cardiac	Heart-LGE [13]	LGE	20	1	Seen
Whole-body	TotalSegmentator-MR [14]	Mixed	263	50	Seen
Hippocampus	MSD-Hippocampus [1]	T1w	260	2	Seen
Vestibular	crossMoDA [15]	T1CE	227	4	Seen
Hip	HipMRI [16]	T2w	211	3	Seen
Spine	SpineMR [17]	T2w	172	15	Seen
Abdominal	AMOS-MR [18]	T1w	60	13	Seen
Shoulder	MSK-Shoulder †	T2w	45	2	Unseen
Gynecology	CervicalCancer [19]	T2w	7	1	Unseen
Total (34)			10,392		26 / 8

G. Per-Dataset and Per-Target Performance

Table S6. DSC by dataset — median (Q1, Q3) across volumes; final row is mean $\pm$ SD.

Domain	Dataset	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2dloss	GVT-v1	GVT-v2	SAMRI-3D
Brain	MSD-BrainTumour	0.24 (0.05, 0.42)	0.26 (0.05, 0.45)	0.23 (0.05, 0.44)	0.21 (0.04, 0.40)	0.23 (0.10, 0.38)	0.34 (0.17, 0.53)	0.57 (0.39, 0.65)	0.33 (0.07, 0.53)	0.40 (0.11, 0.60)	0.31 (0.08, 0.54)	0.40 (0.10, 0.61)
Brain	BraTS-FLAIR	0.82 (0.73, 0.88)	0.81 (0.72, 0.87)	0.84 (0.75, 0.88)	0.83 (0.72, 0.88)	0.71 (0.60, 0.78)	0.76 (0.70, 0.83)	0.84 (0.79, 0.87)	0.89 (0.82, 0.92)	0.89 (0.83, 0.91)	0.90 (0.83, 0.92)	0.89 (0.82, 0.92)
Brain	BraTS-T1	0.71 (0.64, 0.77)	0.72 (0.63, 0.76)	0.71 (0.61, 0.77)	0.68 (0.59, 0.75)	0.59 (0.48, 0.67)	0.60 (0.46, 0.70)	0.76 (0.71, 0.80)	0.79 (0.73, 0.84)	0.79 (0.73, 0.84)	0.81 (0.75, 0.86)	0.80 (0.75, 0.84)
Brain	BraTS-T1CE	0.66 (0.53, 0.80)	0.67 (0.50, 0.80)	0.66 (0.54, 0.79)	0.67 (0.52, 0.80)	0.60 (0.49, 0.76)	0.28 (0.20, 0.40)	0.47 (0.34, 0.62)	0.82 (0.60, 0.90)	0.86 (0.66, 0.91)	0.86 (0.66, 0.91)	0.80 (0.55, 0.89)
Brain	ISLES2022-DWI	0.56 (0.33, 0.83)	0.58 (0.28, 0.75)	0.63 (0.47, 0.78)	0.54 (0.22, 0.78)	0.60 (0.43, 0.73)	0.06 (0.03, 0.18)	0.73 (0.65, 0.78)	0.81 (0.61, 0.86)	0.77 (0.60, 0.86)	0.80 (0.64, 0.88)	0.80 (0.70, 0.88)
Brain	ISLES2022-ADC	0.40 (0.14, 0.59)	0.33 (0.13, 0.60)	0.43 (0.18, 0.59)	0.33 (0.12, 0.61)	0.37 (0.30, 0.52)	0.07 (0.03, 0.15)	0.60 (0.56, 0.67)	0.51 (0.42, 0.71)	0.58 (0.38, 0.72)	0.55 (0.42, 0.73)	0.59 (0.45, 0.71)
Brain	WMH-T1	0.07 (0.03, 0.14)	0.08 (0.04, 0.17)	0.13 (0.06, 0.29)	0.09 (0.04, 0.17)	0.38 (0.25, 0.48)	0.02 (0.01, 0.04)	0.38 (0.29, 0.50)	0.48 (0.32, 0.61)	0.43 (0.25, 0.61)	0.53 (0.40, 0.64)	0.54 (0.40, 0.65)
Brain	WMH-FLAIR	0.13 (0.04, 0.29)	0.11 (0.04, 0.24)	0.17 (0.08, 0.44)	0.15 (0.05, 0.31)	0.42 (0.25, 0.53)	0.02 (0.01, 0.06)	0.45 (0.36, 0.55)	0.59 (0.39, 0.71)	0.46 (0.25, 0.66)	0.58 (0.38, 0.73)	0.54 (0.41, 0.72)
Prostate	PI-CAI	0.81 (0.72, 0.86)	0.78 (0.71, 0.85)	0.78 (0.72, 0.83)	0.77 (0.68, 0.83)	0.70 (0.54, 0.78)	0.24 (0.16, 0.34)	0.88 (0.86, 0.89)	0.92 (0.90, 0.94)	0.91 (0.89, 0.93)	0.92 (0.89, 0.94)	0.91 (0.89, 0.94)
Prostate	Prostate-T2	0.89 (0.86, 0.92)	0.89 (0.85, 0.91)	0.88 (0.85, 0.90)	0.88 (0.85, 0.91)	0.77 (0.60, 0.83)	0.37 (0.28, 0.49)	0.89 (0.88, 0.90)	0.93 (0.91, 0.94)	0.93 (0.91, 0.94)	0.93 (0.90, 0.94)	0.93 (0.91, 0.94)
Prostate	Prostate-ADC	0.83 (0.81, 0.88)	0.85 (0.82, 0.88)	0.85 (0.81, 0.88)	0.86 (0.82, 0.88)	0.78 (0.66, 0.84)	0.76 (0.72, 0.80)	0.89 (0.87, 0.90)	0.93 (0.91, 0.94)	0.93 (0.91, 0.94)	0.92 (0.90, 0.94)	0.92 (0.90, 0.93)
Prostate	QIN-Prostate	0.86 (0.82, 0.89)	0.89 (0.86, 0.90)	0.84 (0.80, 0.90)	0.86 (0.82, 0.90)	0.81 (0.73, 0.86)	0.25 (0.23, 0.53)	0.87 (0.79, 0.89)	0.92 (0.90, 0.93)	0.92 (0.91, 0.93)	0.91 (0.90, 0.92)	0.91 (0.91, 0.92)
Prostate	PROMISE12	0.72 (0.66, 0.76)	0.64 (0.57, 0.71)	0.60 (0.48, 0.70)	0.60 (0.44, 0.65)	0.17 (0.11, 0.21)	0.28 (0.19, 0.37)	0.87 (0.84, 0.89)	0.89 (0.87, 0.90)	0.86 (0.84, 0.88)	0.88 (0.87, 0.89)	0.88 (0.87, 0.90)
Prostate	MSD-Prostate	0.56 (0.18, 0.67)	0.51 (0.20, 0.64)	0.56 (0.24, 0.63)	0.50 (0.21, 0.64)	0.26 (0.22, 0.39)	0.46 (0.24, 0.55)	0.83 (0.50, 0.87)	0.73 (0.34, 0.77)	0.68 (0.28, 0.78)	0.68 (0.26, 0.76)	0.71 (0.34, 0.75)
Prostate	QIN-Lesion	0.65 (0.55, 0.70)	0.62 (0.46, 0.73)	0.54 (0.51, 0.65)	0.62 (0.52, 0.65)	0.60 (0.59, 0.73)	0.35 (0.22, 0.43)	0.62 (0.40, 0.71)	0.70 (0.53, 0.78)	0.68 (0.54, 0.69)	0.65 (0.53, 0.77)	0.66 (0.55, 0.74)
Knee	ZIB-OAI	0.74 (0.55, 0.97)	0.78 (0.53, 0.95)	0.68 (0.45, 0.96)	0.73 (0.48, 0.96)	0.14 (0.01, 0.47)	0.47 (0.01, 0.78)	0.64 (0.33, 0.78)	0.92 (0.79, 0.98)	0.90 (0.79, 0.97)	0.92 (0.80, 0.98)	0.92 (0.78, 0.98)
Knee	OAI-iMorphics-DESS	0.62 (0.45, 0.92)	0.76 (0.48, 0.95)	0.67 (0.37, 0.96)	0.63 (0.39, 0.95)	0.05 (0.02, 0.44)	0.32 (0.15, 0.55)	0.64 (0.48, 0.81)	0.90 (0.80, 0.97)	0.87 (0.79, 0.97)	0.87 (0.80, 0.97)	0.88 (0.78, 0.97)
Knee	OAI-iMorphics-TSE	0.51 (0.32, 0.80)	0.55 (0.37, 0.89)	0.39 (0.22, 0.86)	0.42 (0.26, 0.68)	0.50 (0.17, 0.77)	0.24 (0.10, 0.69)	0.59 (0.43, 0.81)	0.72 (0.57, 0.91)	0.72 (0.56, 0.90)	0.71 (0.58, 0.92)	0.70 (0.54, 0.91)
Knee	OAI-iMorphics-FLASH	0.61 (0.47, 0.90)	0.64 (0.48, 0.93)	0.63 (0.49, 0.94)	0.68 (0.40, 0.92)	0.70 (0.51, 0.88)	0.36 (0.10, 0.50)	0.69 (0.50, 0.80)	0.78 (0.71, 0.97)	0.77 (0.69, 0.94)	0.77 (0.72, 0.95)	0.79 (0.71, 0.96)
Knee	OAI-AKOA	0.61 (0.47, 0.88)	0.71 (0.38, 0.91)	0.54 (0.44, 0.91)	0.61 (0.51, 0.83)	0.46 (0.04, 0.74)	0.47 (0.10, 0.81)	0.73 (0.56, 0.81)	0.88 (0.74, 0.96)	0.89 (0.78, 0.95)	0.88 (0.76, 0.95)	0.87 (0.76, 0.95)
Knee	MSK-T2	0.86 (0.72, 0.94)	0.85 (0.66, 0.93)	0.87 (0.75, 0.94)	0.77 (0.61, 0.93)	0.50 (0.05, 0.73)	0.44 (0.17, 0.82)	0.71 (0.59, 0.82)	0.90 (0.87, 0.96)	0.90 (0.86, 0.95)	0.90 (0.86, 0.95)	0.90 (0.87, 0.95)
Knee	MSK-FLASH	0.18 (0.05, 0.44)	0.12 (0.03, 0.49)	0.11 (0.03, 0.56)	0.08 (0.04, 0.40)	0.55 (0.12, 0.71)	0.40 (0.13, 0.81)	0.74 (0.54, 0.78)	0.38 (0.16, 0.72)	0.50 (0.26, 0.73)	0.68 (0.53, 0.84)	0.71 (0.57, 0.88)
Knee	MSK-PD	0.85 (0.71, 0.92)	0.79 (0.57, 0.91)	0.81 (0.61, 0.92)	0.71 (0.44, 0.90)	0.55 (0.38, 0.74)	0.63 (0.28, 0.83)	0.71 (0.52, 0.79)	0.90 (0.82, 0.95)	0.90 (0.82, 0.95)	0.88 (0.82, 0.94)	0.88 (0.80, 0.94)
Cardiac	ACDC	0.85 (0.79, 0.91)	0.85 (0.79, 0.90)	0.84 (0.78, 0.90)	0.84 (0.76, 0.90)	0.71 (0.49, 0.82)	0.42 (0.28, 0.56)	0.79 (0.57, 0.84)	0.85 (0.80, 0.91)	0.87 (0.82, 0.91)	0.84 (0.78, 0.91)	0.85 (0.80, 0.91)
Cardiac	MSD-Heart	0.50 (0.49, 0.51)	0.44 (0.43, 0.46)	0.44 (0.42, 0.46)	0.41 (0.39, 0.43)	0.83 (0.82, 0.84)	0.71 (0.71, 0.72)	0.80 (0.79, 0.81)	0.87 (0.86, 0.87)	0.85 (0.84, 0.86)	0.87 (0.87, 0.87)	0.89 (0.88, 0.89)
Cardiac	Heart-LGE	0.69 (0.66, 0.72)	0.74 (0.73, 0.76)	0.68 (0.66, 0.70)	0.72 (0.72, 0.72)	0.74 (0.69, 0.78)	0.71 (0.71, 0.72)	0.82 (0.82, 0.82)	0.87 (0.85, 0.88)	0.88 (0.87, 0.89)	0.87 (0.86, 0.89)	0.87 (0.86, 0.89)
Whole-body	TotalSegmentator-MR	0.23 (0.06, 0.52)	0.23 (0.05, 0.55)	0.17 (0.02, 0.49)	0.17 (0.02, 0.45)	0.36 (0.17, 0.59)	0.24 (0.10, 0.46)	0.66 (0.48, 0.77)	0.52 (0.26, 0.76)	0.47 (0.15, 0.75)	0.55 (0.28, 0.77)	0.57 (0.31, 0.77)
Hippocampus	MSD-Hippocampus	0.66 (0.55, 0.74)	0.65 (0.51, 0.75)	0.58 (0.50, 0.66)	0.63 (0.50, 0.70)	0.05 (0.02, 0.10)	0.20 (0.16, 0.24)	0.74 (0.67, 0.76)	0.85 (0.82, 0.87)	0.85 (0.80, 0.88)	0.84 (0.80, 0.87)	0.84 (0.80, 0.86)
Vestibular	crossMoDA	0.37 (0.22, 0.64)	0.34 (0.20, 0.57)	0.33 (0.15, 0.66)	0.42 (0.19, 0.65)	0.47 (0.35, 0.62)	0.22 (0.09, 0.43)	0.45 (0.25, 0.71)	0.60 (0.53, 0.80)	0.62 (0.46, 0.82)	0.63 (0.51, 0.81)	0.63 (0.49, 0.79)
Hip	HipMRI	0.64 (0.39, 0.87)	0.55 (0.26, 0.75)	0.61 (0.42, 0.78)	0.48 (0.28, 0.69)	0.48 (0.37, 0.67)	0.58 (0.33, 0.73)	0.74 (0.72, 0.78)	0.79 (0.71, 0.95)	0.77 (0.66, 0.95)	0.83 (0.75, 0.95)	0.87 (0.77, 0.95)
Spine	SpineMR	0.81 (0.74, 0.84)	0.78 (0.71, 0.82)	0.80 (0.74, 0.84)	0.80 (0.73, 0.84)	0.48 (0.18, 0.72)	0.28 (0.09, 0.55)	0.65 (0.60, 0.69)	0.86 (0.83, 0.88)	0.86 (0.83, 0.88)	0.86 (0.83, 0.88)	0.84 (0.82, 0.87)
Abdominal	AMOS-MR	0.38 (0.08, 0.69)	0.28 (0.07, 0.64)	0.35 (0.07, 0.72)	0.24 (0.07, 0.72)	0.40 (0.27, 0.64)	0.28 (0.13, 0.49)	0.72 (0.61, 0.80)	0.71 (0.42, 0.89)	0.64 (0.35, 0.88)	0.72 (0.47, 0.88)	0.70 (0.47, 0.87)
Shoulder	MSK-Shoulder	0.36 (0.00, 0.53)	0.36 (0.07, 0.48)	0.17 (0.00, 0.40)	0.36 (0.07, 0.46)	0.24 (0.02, 0.39)	0.14 (0.04, 0.28)	0.27 (0.14, 0.56)	0.59 (0.12, 0.75)	0.63 (0.56, 0.74)	0.44 (0.06, 0.76)	0.72 (0.07, 0.77)
Gynecology	CervicalCancer	0.70 (0.70, 0.70)	0.67 (0.67, 0.67)	0.63 (0.63, 0.63)	0.66 (0.66, 0.66)	0.53 (0.53, 0.53)	0.55 (0.55, 0.55)	0.68 (0.68, 0.68)	0.74 (0.74, 0.74)	0.75 (0.75, 0.75)	0.73 (0.73, 0.73)	0.76 (0.76, 0.76)
Mean $\pm$ SD		0.57 $\pm$ 0.23	0.58 $\pm$ 0.24	0.56 $\pm$ 0.24	0.56 $\pm$ 0.24	0.49 $\pm$ 0.22	0.37 $\pm$ 0.21	0.69 $\pm$ 0.15	0.76 $\pm$ 0.17	0.76 $\pm$ 0.16	0.77 $\pm$ 0.16	0.78 $\pm$ 0.14

Table S7. HD by dataset — median (Q1, Q3) across volumes; final row is mean $\pm$ SD.

Domain	Dataset	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2dloss	GVT-v1	GVT-v2	SAMRI-3D
Brain	MSD-BrainTumour	63.90 (53.55, 76.05)	59.77 (43.15, 72.10)	57.22 (41.21, 71.12)	65.72 (50.28, 77.99)	48.43 (32.16, 62.45)	21.47 (15.02, 29.70)	14.46 (9.26, 21.16)	35.28 (21.22, 49.45)	26.48 (17.81, 41.86)	40.56 (23.40, 64.26)	38.41 (21.62, 52.73)
Brain	BraTS-FLAIR	22.19 (14.04, 30.95)	24.51 (15.69, 37.59)	17.84 (11.79, 34.77)	24.93 (16.06, 40.72)	30.13 (20.17, 39.99)	13.96 (10.92, 19.40)	19.44 (11.57, 53.37)	14.05 (8.77, 24.62)	14.92 (9.68, 22.60)	13.25 (8.46, 25.81)	13.38 (8.53, 21.85)
Brain	BraTS-T1	24.71 (17.51, 35.07)	29.91 (19.80, 37.10)	29.05 (20.47, 38.46)	34.43 (26.21, 51.98)	27.08 (23.33, 34.02)	20.00 (15.57, 26.21)	17.01 (11.97, 21.45)	18.07 (13.42, 26.81)	18.03 (14.40, 27.24)	19.40 (12.61, 31.76)	16.76 (12.78, 23.19)
Brain	BraTS-T1CE	8.80 (6.32, 11.58)	9.27 (6.65, 13.38)	8.80 (6.71, 12.44)	8.63 (6.10, 12.28)	9.77 (7.35, 16.26)	11.18 (9.43, 15.20)	8.15 (6.16, 12.53)	7.62 (5.41, 10.68)	7.07 (5.20, 10.63)	7.45 (5.02, 10.04)	7.75 (5.54, 11.52)
Brain	ISLES2022-DWI	15.78 (8.31, 26.40)	17.42 (10.46, 25.51)	11.49 (4.36, 24.10)	19.10 (9.50, 26.94)	5.61 (2.24, 9.84)	16.03 (13.03, 19.12)	2.83 (2.00, 4.21)	3.74 (1.49, 7.45)	3.87 (2.06, 7.88)	3.24 (1.41, 5.17)	3.32 (1.73, 5.92)
Brain	ISLES2022-ADC	16.58 (5.10, 26.60)	20.55 (8.00, 24.57)	12.08 (4.86, 18.69)	22.13 (8.97, 29.00)	10.03 (5.49, 17.22)	16.36 (14.01, 19.67)	3.16 (2.40, 5.02)	4.69 (2.93, 8.65)	4.58 (3.12, 8.80)	4.47 (2.96, 7.98)	4.24 (2.45, 7.87)
Brain	WMH-T1	45.84 (32.75, 62.95)	44.37 (32.83, 66.59)	32.77 (14.46, 53.24)	50.98 (33.83, 68.88)	6.70 (3.32, 16.44)	16.42 (14.06, 20.55)	5.00 (3.61, 8.99)	6.78 (3.46, 11.26)	8.60 (3.94, 14.11)	4.53 (2.40, 8.08)	4.58 (2.45, 9.26)
Brain	WMH-FLAIR	55.18 (33.77, 83.63)	54.84 (35.56, 82.97)	48.21 (23.17, 74.44)	56.87 (37.43, 84.24)	8.18 (3.74, 20.48)	16.82 (14.18, 21.03)	6.00 (3.61, 44.29)	7.78 (4.12, 17.67)	15.23 (5.81, 27.65)	6.56 (3.04, 16.75)	6.00 (3.32, 16.20)
Prostate	PI-CAI	34.07 (15.61, 53.84)	40.14 (24.41, 60.61)	41.78 (22.73, 69.01)	56.59 (34.00, 87.18)	45.11 (30.19, 60.59)	40.52 (34.21, 50.68)	14.63 (11.67, 17.66)	6.56 (4.69, 8.86)	8.15 (5.83, 11.69)	8.03 (5.05, 13.81)	8.22 (5.17, 13.18)
Prostate	Prostate-T2	11.38 (7.26, 15.51)	11.46 (7.03, 18.71)	13.15 (8.11, 17.90)	14.88 (8.87, 21.41)	23.16 (13.71, 40.95)	35.90 (25.53, 44.73)	13.15 (10.51, 15.30)	6.00 (5.02, 7.12)	6.04 (5.02, 7.14)	6.12 (5.12, 8.02)	5.92 (5.00, 7.07)
Prostate	Prostate-ADC	9.38 (5.10, 12.33)	6.16 (5.10, 12.33)	8.60 (4.58, 10.82)	7.28 (4.47, 12.08)	8.72 (5.92, 12.53)	7.21 (6.16, 8.60)	2.83 (2.45, 3.61)	2.45 (2.24, 3.00)	2.83 (2.24, 3.46)	2.83 (2.24, 3.61)	2.83 (2.24, 3.61)
Prostate	QIN-Prostate	18.06 (10.00, 32.16)	11.18 (9.49, 18.79)	24.54 (10.25, 42.56)	13.34 (11.75, 33.66)	24.56 (18.14, 29.53)	37.54 (22.16, 70.72)	11.58 (10.63, 14.32)	7.14 (7.00, 9.90)	9.49 (6.40, 10.49)	8.60 (7.68, 20.10)	7.87 (6.48, 9.49)
Prostate	PROMISE12	47.98 (42.36, 87.95)	50.79 (48.72, 91.99)	86.98 (56.81, 116.98)	86.27 (69.00, 100.63)	35.83 (11.54, 43.55)	25.83 (22.24, 54.64)	10.77 (6.52, 14.23)	6.48 (4.74, 9.72)	10.00 (8.71, 17.49)	5.92 (4.79, 8.20)	7.07 (6.08, 9.35)
Prostate	MSD-Prostate	49.17 (29.61, 56.39)	41.56 (29.70, 52.54)	45.61 (23.34, 58.98)	39.10 (23.11, 46.11)	42.59 (33.67, 51.02)	20.05 (16.61, 22.63)	9.87 (8.11, 13.04)	8.68 (7.28, 15.32)	11.38 (7.59, 13.33)	9.72 (7.93, 17.50)	10.74 (7.35, 16.18)
Prostate	QIN-Lesion	9.43 (5.73, 17.18)	9.00 (4.03, 12.04)	9.54 (6.53, 12.26)	7.35 (4.87, 10.46)	8.60 (4.96, 11.52)	12.57 (11.52, 13.74)	7.87 (5.54, 10.00)	5.92 (3.86, 10.23)	5.92 (3.99, 10.14)	5.74 (4.30, 9.51)	7.87 (2.99, 11.25)
Knee/MSK	ZIB-OAI	83.93 (40.03, 141.11)	82.89 (57.93, 130.45)	77.81 (55.33, 116.48)	104.78 (69.05, 141.56)	78.72 (62.76, 92.44)	46.50 (37.96, 76.59)	52.98 (30.83, 70.26)	14.95 (10.59, 24.14)	18.57 (12.44, 40.29)	14.51 (10.42, 22.13)	20.40 (13.98, 33.64)
Knee/MSK	OAI-iMorphics-DESS	127.80 (59.45, 203.48)	99.56 (47.78, 166.77)	76.30 (30.96, 123.98)	118.88 (80.14, 170.69)	68.21 (53.98, 86.64)	57.16 (39.46, 70.40)	39.94 (16.48, 61.63)	13.34 (8.22, 21.31)	16.14 (9.84, 26.25)	13.42 (8.69, 20.69)	19.03 (9.70, 28.27)
Knee/MSK	OAI-iMorphics-TSE	74.77 (41.02, 108.06)	70.47 (40.02, 105.53)	69.57 (45.45, 101.33)	112.75 (76.64, 170.42)	51.10 (36.68, 68.47)	48.76 (37.54, 56.68)	38.33 (15.24, 56.32)	21.66 (14.43, 39.88)	23.97 (16.05, 38.02)	21.67 (14.80, 35.46)	23.02 (15.73, 44.73)
Knee/MSK	OAI-iMorphics-FLASH	176.21 (117.75, 245.24)	103.42 (59.08, 187.78)	63.65 (24.72, 155.07)	145.42 (88.86, 221.05)	61.06 (37.46, 90.59)	52.48 (37.31, 69.63)	58.21 (31.54, 75.64)	13.94 (11.12, 26.14)	19.70 (12.40, 51.47)	22.11 (11.51, 36.72)	18.19 (11.12, 46.24)
Knee/MSK	OAI-AKOA	94.45 (60.04, 174.93)	88.88 (60.62, 148.13)	67.71 (58.37, 97.94)	108.02 (76.02, 199.64)	59.21 (49.73, 79.23)	42.96 (32.14, 69.24)	27.22 (12.18, 48.85)	12.80 (8.02, 39.42)	13.82 (9.49, 20.91)	8.74 (7.33, 15.07)	18.13 (9.24, 36.97)
Knee/MSK	MSK-T2	21.00 (7.93, 80.33)	65.51 (16.04, 97.46)	15.56 (9.75, 56.11)	91.68 (57.47, 147.79)	63.70 (33.11, 94.88)	35.73 (25.89, 46.22)	33.43 (12.96, 60.32)	10.86 (8.48, 40.34)	13.19 (8.72, 46.89)	17.30 (8.40, 58.40)	21.79 (8.88, 72.68)
Knee/MSK	MSK-FLASH	104.79 (64.08, 156.99)	140.28 (68.96, 156.90)	111.32 (46.60, 175.74)	136.71 (80.25, 151.90)	52.10 (37.05, 121.62)	33.04 (28.92, 54.86)	30.35 (16.44, 58.72)	52.94 (33.52, 116.23)	40.92 (16.61, 104.42)	49.48 (30.52, 115.99)	35.37 (18.01, 129.30)
Knee/MSK	MSK-PD	51.62 (17.96, 80.09)	48.96 (30.72, 78.69)	43.34 (14.37, 78.81)	94.61 (62.53, 123.78)	37.45 (17.30, 64.96)	45.46 (25.18, 54.02)	42.26 (24.29, 73.04)	12.22 (10.56, 18.29)	12.77 (10.05, 74.37)	14.12 (11.04, 19.99)	19.80 (11.38, 48.02)
Cardiac	ACDC	6.67 (4.41, 10.05)	6.75 (4.44, 10.55)	7.07 (4.90, 11.36)	8.06 (4.92, 11.88)	8.63 (6.08, 12.75)	33.13 (17.83, 52.32)	8.54 (5.79, 11.22)	4.79 (3.00, 7.07)	4.58 (3.00, 6.82)	5.00 (3.24, 8.08)	4.90 (3.16, 7.02)
Cardiac	MSD-Heart	67.12 (65.75, 68.48)	70.77 (70.16, 71.39)	67.82 (66.73, 68.91)	70.52 (69.76, 71.28)	15.18 (14.70, 15.67)	22.55 (18.71, 26.39)	10.72 (10.56, 10.88)	17.79 (15.00, 20.58)	28.02 (24.27, 31.77)	21.06 (19.26, 22.86)	13.33 (12.77, 13.89)
Cardiac	Heart-LGE	49.76 (43.57, 55.95)	41.44 (39.17, 43.71)	38.95 (37.93, 39.98)	36.56 (36.39, 36.73)	23.26 (21.57, 24.95)	17.07 (15.83, 18.31)	73.58 (73.32, 73.85)	14.26 (12.15, 16.36)	18.75 (18.61, 18.89)	13.14 (11.89, 14.40)	12.19 (11.41, 12.97)
Whole-body	TotalSegmentator-MR	78.42 (39.38, 122.70)	70.18 (40.83, 108.69)	73.82 (39.89, 116.57)	94.37 (56.67, 132.36)	41.74 (22.42, 67.87)	30.35 (18.14, 49.42)	27.62 (10.12, 76.12)	26.44 (14.28, 47.97)	30.84 (15.17, 56.18)	28.81 (13.75, 53.82)	27.69 (13.92, 49.96)
Hippocampus	MSD-Hippocampus	13.64 (9.19, 18.98)	13.75 (9.35, 19.54)	15.20 (12.81, 19.90)	15.65 (12.88, 21.31)	12.86 (10.05, 14.47)	16.20 (13.96, 19.53)	4.24 (3.57, 7.11)	3.67 (3.00, 4.24)	4.06 (3.16, 5.02)	3.61 (3.00, 4.24)	3.61 (3.12, 4.24)
Vestibular	crossMoDA	15.23 (6.65, 27.51)	19.84 (9.11, 39.52)	18.91 (11.72, 40.43)	19.71 (8.77, 32.36)	5.83 (4.18, 8.76)	9.67 (7.63, 11.49)	32.11 (7.09, 99.02)	4.57 (3.16, 6.40)	5.39 (3.16, 7.66)	4.18 (3.00, 5.72)	4.24 (3.00, 6.16)
Hip	HipMRI	35.93 (10.89, 57.69)	51.02 (25.21, 79.15)	51.81 (32.77, 75.47)	81.41 (66.77, 91.61)	57.64 (39.16, 75.58)	17.07 (11.72, 23.36)	5.83 (4.69, 10.40)	20.28 (8.56, 39.47)	23.95 (8.38, 40.28)	15.54 (8.75, 30.56)	13.69 (7.31, 29.38)
Spine	SpineMR	28.07 (10.95, 53.39)	37.12 (13.60, 64.15)	34.55 (11.40, 68.45)	38.20 (12.17, 75.34)	53.08 (16.16, 78.32)	57.76 (21.66, 99.03)	33.20 (23.11, 53.11)	11.66 (6.32, 20.90)	12.04 (6.40, 23.11)	11.18 (6.78, 22.02)	14.76 (7.00, 30.55)
Abdominal	AMOS-MR	58.69 (28.51, 76.83)	57.79 (33.97, 96.22)	56.41 (30.28, 91.54)	68.18 (45.17, 109.59)	37.68 (21.96, 52.04)	20.83 (12.67, 33.53)	14.69 (7.83, 83.86)	14.46 (6.54, 29.64)	17.76 (8.70, 39.59)	14.85 (6.54, 35.77)	15.23 (6.42, 31.98)
Shoulder	MSK-Shoulder	107.85 (63.69, 126.07)	128.54 (112.04, 185.90)	117.62 (96.41, 160.40)	156.58 (134.98, 192.33)	60.13 (49.64, 75.93)	75.93 (53.71, 95.87)	62.61 (20.32, 86.72)	54.95 (26.14, 114.18)	57.10 (28.01, 98.38)	45.29 (13.89, 99.49)	23.80 (14.87, 97.86)
Gynecology	CervicalCancer	17.61 (17.61, 17.61)	15.81 (15.81, 15.81)	12.53 (12.53, 12.53)	16.25 (16.25, 16.25)	17.12 (17.12, 17.12)	13.04 (13.04, 13.04)	14.87 (14.87, 14.87)	16.19 (16.19, 16.19)	15.84 (15.84, 15.84)	15.84 (15.84, 15.84)	16.79 (16.79, 16.79)
Mean $\pm$ SD		48.41 $\pm$ 39.64	48.35 $\pm$ 35.14	43.16 $\pm$ 30.30	59.59 $\pm$ 44.37	33.51 $\pm$ 21.81	29.05 $\pm$ 16.70	22.28 $\pm$ 18.82	14.21 $\pm$ 12.32	15.59 $\pm$ 11.43	14.30 $\pm$ 11.68	13.85 $\pm$ 9.01

Table S8. MSD by dataset — median (Q1, Q3) across volumes; final row is mean $\pm$ SD.

Domain	Dataset	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2dloss	GVT-v1	GVT-v2	SAMRI-3D
Brain	MSD-BrainTumour	12.26 (8.70, 16.41)	10.54 (6.99, 14.59)	9.89 (6.74, 14.34)	11.50 (7.99, 16.33)	8.85 (6.02, 12.52)	4.08 (3.02, 5.87)	2.48 (2.04, 3.21)	6.39 (4.00, 9.04)	4.92 (3.22, 7.06)	6.30 (3.94, 9.23)	5.67 (3.55, 8.27)
Brain	BraTS-FLAIR	3.32 (2.54, 4.55)	3.50 (2.58, 4.97)	3.18 (2.23, 4.36)	3.26 (2.20, 4.94)	4.67 (3.48, 5.96)	2.36 (2.05, 2.79)	2.42 (2.08, 2.91)	2.24 (1.74, 2.94)	2.24 (1.70, 3.02)	1.93 (1.60, 2.64)	2.07 (1.65, 3.07)
Brain	BraTS-T1	4.45 (3.52, 5.82)	4.64 (3.73, 6.01)	4.83 (3.60, 6.04)	5.23 (3.81, 7.34)	5.50 (4.44, 7.11)	4.13 (3.57, 5.23)	3.02 (2.60, 3.45)	3.29 (2.60, 4.29)	3.46 (2.58, 4.32)	2.88 (2.46, 3.75)	3.31 (2.49, 4.02)
Brain	BraTS-T1CE	2.32 (1.67, 3.18)	2.62 (1.68, 3.27)	2.31 (1.69, 3.31)	2.37 (1.78, 3.29)	2.48 (1.94, 3.13)	3.14 (2.70, 3.73)	1.90 (1.71, 2.22)	1.68 (1.32, 2.46)	1.55 (1.26, 2.15)	1.47 (1.25, 2.40)	1.67 (1.32, 2.74)
Brain	ISLES2022-DW1	3.72 (1.61, 6.17)	3.56 (2.20, 6.30)	3.12 (1.37, 5.87)	4.95 (1.86, 6.41)	1.51 (1.13, 2.98)	6.25 (4.78, 7.63)	1.16 (1.10, 1.26)	1.29 (1.04, 1.82)	1.32 (1.07, 1.82)	1.17 (1.03, 1.51)	1.27 (1.07, 1.59)
Brain	ISLES2022-ADC	3.25 (1.49, 5.49)	4.26 (1.97, 6.50)	2.88 (1.57, 4.23)	3.85 (2.05, 6.62)	2.82 (1.80, 3.64)	5.28 (4.13, 6.66)	1.22 (1.16, 1.40)	1.49 (1.17, 2.33)	1.60 (1.18, 2.22)	1.49 (1.19, 2.45)	1.41 (1.17, 2.21)
Brain	WMH-T1	9.01 (5.90, 12.37)	8.30 (5.70, 12.39)	6.72 (3.45, 11.17)	9.91 (6.22, 15.05)	2.19 (1.32, 4.84)	6.90 (5.94, 7.81)	1.55 (1.33, 2.09)	1.85 (1.30, 2.62)	2.04 (1.34, 3.26)	1.48 (1.18, 2.08)	1.49 (1.22, 2.30)
Brain	WMH-FLAIR	9.79 (5.72, 17.75)	8.85 (5.37, 13.06)	7.09 (3.41, 14.01)	10.20 (6.08, 16.12)	2.21 (1.30, 4.55)	6.84 (5.70, 7.94)	1.67 (1.35, 2.08)	1.67 (1.23, 3.39)	2.77 (1.43, 5.11)	1.63 (1.18, 2.92)	1.57 (1.19, 2.79)
Prostate	PI-CAI	3.42 (2.36, 4.74)	3.79 (3.09, 5.32)	4.08 (3.09, 6.00)	4.84 (3.75, 7.74)	4.61 (3.28, 5.70)	9.24 (7.78, 11.02)	2.04 (1.74, 2.19)	1.49 (1.29, 1.64)	1.68 (1.40, 2.01)	1.58 (1.46, 1.89)	1.58 (1.41, 1.79)
Prostate	Prostate-T2	2.09 (1.75, 2.50)	2.43 (1.93, 2.94)	2.39 (1.91, 2.75)	2.37 (2.12, 2.88)	3.35 (2.65, 5.06)	9.19 (7.16, 11.87)	2.00 (1.79, 2.19)	1.47 (1.38, 1.73)	1.48 (1.33, 1.86)	1.48 (1.38, 1.75)	1.50 (1.38, 1.62)
Prostate	Prostate-ADC	1.66 (1.39, 1.95)	1.56 (1.32, 1.83)	1.58 (1.42, 1.75)	1.50 (1.32, 1.72)	1.94 (1.56, 2.52)	2.21 (1.96, 2.42)	1.29 (1.19, 1.55)	1.12 (1.08, 1.28)	1.14 (1.07, 1.25)	1.15 (1.08, 1.32)	1.14 (1.10, 1.25)
Prostate	QIN-Prostate	3.42 (3.01, 3.66)	2.51 (2.04, 3.57)	3.86 (2.85, 5.11)	3.02 (2.53, 3.66)	2.98 (2.57, 4.21)	13.08 (5.46, 17.60)	2.37 (2.29, 3.86)	1.58 (1.53, 2.50)	1.63 (1.57, 2.75)	1.91 (1.59, 2.80)	1.64 (1.51, 2.43)
Prostate	PROMISE12	4.97 (3.46, 7.85)	5.82 (4.96, 8.29)	8.83 (6.66, 12.46)	8.22 (7.07, 10.64)	4.77 (3.29, 6.78)	6.55 (5.49, 10.95)	1.72 (1.58, 1.88)	1.45 (1.34, 1.52)	1.80 (1.58, 2.19)	1.45 (1.33, 1.66)	1.47 (1.40, 1.78)
Prostate	MSD-Prostate	5.57 (4.23, 5.97)	5.12 (4.02, 6.19)	5.68 (3.51, 7.52)	5.05 (3.87, 5.78)	4.95 (4.31, 7.39)	4.48 (3.62, 5.38)	2.42 (1.63, 2.68)	2.12 (1.85, 2.79)	2.35 (2.02, 2.99)	2.30 (2.17, 2.67)	2.23 (1.93, 2.49)
Prostate	QIN-Lesion	2.13 (1.43, 3.19)	2.17 (1.41, 2.44)	2.63 (2.03, 2.99)	2.11 (1.62, 2.14)	2.24 (1.56, 2.47)	6.14 (4.90, 6.60)	2.36 (1.44, 3.18)	1.51 (1.29, 1.90)	1.53 (1.41, 1.84)	1.51 (1.37, 1.83)	1.74 (1.34, 2.19)
Knee/MSK	ZIB-OAI	8.35 (3.48, 13.03)	8.09 (5.11, 11.79)	8.94 (5.01, 12.90)	9.80 (5.50, 14.25)	19.25 (13.66, 24.61)	8.98 (6.91, 15.06)	7.51 (5.65, 9.58)	1.71 (1.43, 2.56)	1.97 (1.51, 3.24)	1.66 (1.42, 2.39)	2.03 (1.57, 2.91)
Knee/MSK	OAI-iMorphics-DESS	12.43 (4.30, 35.59)	9.32 (3.81, 22.23)	8.04 (2.69, 14.78)	10.05 (5.35, 20.34)	16.15 (10.62, 24.13)	11.53 (7.33, 16.27)	4.95 (2.86, 7.66)	1.60 (1.34, 2.30)	1.68 (1.35, 2.84)	1.56 (1.32, 2.13)	1.81 (1.42, 3.02)
Knee/MSK	OAI-iMorphics-TSE	9.83 (4.34, 15.66)	7.16 (4.07, 12.19)	8.44 (5.58, 14.47)	14.11 (8.44, 23.04)	9.66 (5.92, 15.10)	8.61 (6.85, 11.09)	4.54 (2.90, 6.91)	2.84 (2.17, 5.30)	2.92 (2.33, 5.04)	2.80 (2.13, 4.71)	3.05 (2.28, 5.90)
Knee/MSK	OAI-iMorphics-FLASH	19.03 (11.20, 41.19)	7.98 (5.35, 22.37)	7.08 (1.85, 15.00)	12.21 (7.28, 37.54)	7.97 (3.95, 16.65)	9.34 (5.95, 15.67)	3.83 (2.92, 7.02)	1.80 (1.43, 2.45)	1.88 (1.53, 4.65)	2.00 (1.59, 3.30)	2.06 (1.56, 4.47)
Knee/MSK	OAI-AKOA	9.49 (3.94, 34.13)	10.48 (7.63, 25.50)	11.23 (6.10, 17.19)	10.89 (8.28, 17.51)	12.09 (9.21, 15.39)	8.28 (5.60, 13.21)	4.28 (2.64, 6.09)	1.59 (1.32, 4.28)	1.59 (1.34, 2.96)	1.43 (1.33, 1.92)	1.57 (1.42, 4.21)
Knee/MSK	MSK-T2	2.81 (1.50, 8.23)	4.91 (2.27, 10.78)	2.36 (1.84, 4.21)	9.76 (4.11, 19.48)	10.65 (7.04, 16.59)	6.38 (4.51, 8.76)	3.98 (2.53, 7.24)	1.67 (1.31, 3.54)	1.76 (1.40, 5.87)	2.01 (1.32, 5.62)	2.29 (1.41, 8.44)
Knee/MSK	MSK-FLASH	18.67 (9.99, 31.51)	23.09 (8.71, 28.24)	13.52 (6.34, 31.52)	23.38 (11.55, 36.07)	8.96 (5.71, 19.13)	7.23 (5.20, 9.75)	4.33 (3.27, 6.51)	7.34 (3.81, 20.06)	5.09 (2.15, 19.87)	4.59 (2.98, 14.32)	4.29 (2.41, 7.85)
Knee/MSK	MSK-PD	2.49 (2.19, 8.25)	4.22 (2.52, 9.83)	3.02 (1.68, 8.01)	8.54 (4.81, 14.96)	6.45 (2.59, 10.72)	5.50 (4.33, 7.72)	4.27 (2.79, 6.71)	1.79 (1.40, 1.99)	1.77 (1.43, 3.38)	1.87 (1.58, 2.09)	1.93 (1.70, 3.64)
Cardiac	ACDC	1.51 (1.26, 1.86)	1.50 (1.27, 1.85)	1.54 (1.24, 1.94)	1.59 (1.27, 2.05)	1.91 (1.56, 2.48)	4.49 (3.32, 5.88)	1.74 (1.50, 2.06)	1.37 (1.14, 1.68)	1.35 (1.12, 1.64)	1.38 (1.17, 1.69)	1.33 (1.14, 1.63)
Cardiac	MSD-Heart	11.47 (11.01, 11.93)	12.35 (12.21, 12.49)	11.64 (11.38, 11.90)	12.85 (12.82, 12.88)	2.57 (2.48, 2.67)	4.02 (3.78, 4.26)	2.03 (2.01, 2.04)	2.66 (2.35, 2.97)	3.10 (2.74, 3.45)	2.85 (2.57, 3.13)	2.12 (2.08, 2.15)
Cardiac	Heart-LGE	6.45 (5.29, 7.60)	3.91 (3.50, 4.32)	5.66 (5.01, 6.31)	4.31 (4.24, 4.38)	3.77 (3.46, 4.09)	3.77 (3.59, 3.95)	2.09 (2.03, 2.16)	1.93 (1.74, 2.12)	1.75 (1.74, 1.76)	1.90 (1.76, 2.04)	1.93 (1.79, 2.06)
Whole-body	TotalSegmentator-MR	11.63 (5.08, 22.48)	10.10 (5.21, 18.45)	11.53 (5.03, 21.75)	14.54 (6.41, 24.45)	6.93 (3.80, 12.80)	6.68 (4.28, 11.03)	3.03 (1.91, 5.17)	4.09 (2.44, 7.85)	4.58 (2.58, 9.18)	4.08 (2.29, 7.86)	3.97 (2.38, 7.78)
Hippocampus	MSD-Hippocampus	2.32 (1.77, 3.14)	2.54 (1.86, 3.60)	2.77 (2.07, 3.46)	2.72 (2.22, 3.67)	3.39 (2.93, 3.91)	5.06 (4.42, 5.73)	1.36 (1.28, 1.70)	1.19 (1.13, 1.29)	1.21 (1.15, 1.36)	1.21 (1.15, 1.30)	1.21 (1.15, 1.29)
Vestibular	crossMoDA	2.80 (1.66, 5.62)	3.59 (2.27, 6.12)	3.60 (2.37, 5.91)	3.05 (1.79, 5.30)	1.71 (1.41, 2.18)	3.49 (2.76, 5.25)	2.51 (1.63, 4.15)	1.44 (1.20, 1.80)	1.43 (1.19, 2.04)	1.36 (1.18, 1.56)	1.37 (1.23, 1.62)
Hip	HipMRI	6.27 (1.96, 10.47)	7.91 (3.69, 11.88)	7.26 (4.37, 13.77)	13.95 (8.63, 19.16)	9.36 (5.68, 15.29)	4.33 (2.74, 5.57)	1.89 (1.57, 2.36)	2.82 (1.41, 6.52)	4.06 (1.46, 7.52)	2.45 (1.40, 4.61)	1.90 (1.33, 3.91)
Spine	SpineMR	2.80 (1.80, 5.01)	3.22 (2.11, 5.56)	3.22 (1.91, 6.86)	3.49 (1.93, 7.30)	6.62 (2.65, 14.58)	9.05 (3.91, 20.09)	4.23 (3.49, 5.14)	1.80 (1.45, 2.24)	1.76 (1.46, 2.29)	1.75 (1.44, 2.23)	1.93 (1.55, 2.61)
Abdominal	AMOS-MR	9.00 (3.80, 14.67)	8.35 (5.14, 12.22)	9.87 (4.18, 16.43)	10.45 (5.43, 17.10)	6.99 (3.83, 11.31)	5.16 (3.70, 8.01)	2.32 (1.69, 3.28)	2.51 (1.49, 4.63)	2.67 (1.78, 5.69)	2.41 (1.58, 5.06)	2.40 (1.54, 4.69)
Shoulder	MSK-Shoulder	17.83 (9.25, 19.39)	14.08 (11.76, 23.86)	20.48 (13.71, 32.28)	19.72 (12.84, 31.73)	9.79 (5.74, 15.27)	12.18 (6.49, 15.54)	7.93 (2.90, 11.26)	8.30 (2.76, 14.36)	4.99 (3.26, 10.32)	7.35 (2.48, 17.53)	2.74 (2.42, 18.11)
Gynecology	CervicalCancer	2.29 (2.29, 2.29)	2.39 (2.39, 2.39)	2.57 (2.57, 2.57)	2.40 (2.40, 2.40)	1.99 (1.99, 1.99)	3.28 (3.28, 3.28)	2.47 (2.47, 2.47)	2.12 (2.12, 2.12)	2.08 (2.08, 2.08)	1.98 (1.98, 1.98)	2.06 (2.06, 2.06)
Mean $\pm$ SD		6.73 $\pm$ 5.07	6.32 $\pm$ 4.46	6.23 $\pm$ 4.22	7.83 $\pm$ 5.44	5.92 $\pm$ 4.26	6.39 $\pm$ 2.78	2.91 $\pm$ 1.62	2.39 $\pm$ 1.70	2.33 $\pm$ 1.15	2.25 $\pm$ 1.39	2.11 $\pm$ 0.96

Table S9. DSC by target — median (Q1, Q3) across volumes; final row is mean $\pm$ SD.

Dataset	Target	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMR1-3D-2dloss	GVT-v1	GVT-v2	SAMR1-3D
MSD-BrainTumour	Edema	0.31 (0.17,0.49)	0.33 (0.19,0.49)	0.30 (0.15,0.46)	0.29 (0.16,0.45)	0.24 (0.12,0.38)	0.43 (0.27,0.62)	0.57 (0.44,0.63)	0.35 (0.19,0.53)	0.44 (0.27,0.61)	0.39 (0.22,0.57)	0.45 (0.28,0.62)
MSD-BrainTumour	Enhancing tumour	0.35 (0.23,0.52)	0.40 (0.27,0.55)	0.42 (0.24,0.59)	0.35 (0.22,0.51)	0.25 (0.11,0.39)	0.39 (0.17,0.55)	0.63 (0.55,0.69)	0.59 (0.36,0.63)	0.57 (0.44,0.68)	0.52 (0.36,0.62)	0.56 (0.43,0.69)
MSD-BrainTumour	Non-enhancing tumour	0.02 (0.00,0.08)	0.03 (0.00,0.08)	0.02 (0.00,0.09)	0.02 (0.00,0.06)	0.19 (0.08,0.36)	0.22 (0.12,0.37)	0.44 (0.30,0.61)	0.03 (0.00,0.11)	0.05 (0.00,0.12)	0.04 (0.01,0.11)	0.06 (0.01,0.18)
BraTS-FLAIR	Brain Tumor	0.82 (0.73,0.88)	0.81 (0.72,0.87)	0.84 (0.75,0.88)	0.83 (0.72,0.88)	0.71 (0.60,0.78)	0.76 (0.70,0.83)	0.84 (0.79,0.87)	0.89 (0.82,0.92)	0.90 (0.83,0.91)	0.90 (0.83,0.92)	0.89 (0.82,0.92)
BraTS-T1	Brain Tumor	0.71 (0.64,0.77)	0.72 (0.63,0.76)	0.71 (0.61,0.77)	0.68 (0.59,0.75)	0.59 (0.48,0.67)	0.60 (0.46,0.70)	0.76 (0.71,0.80)	0.79 (0.73,0.84)	0.79 (0.73,0.84)	0.81 (0.75,0.86)	0.80 (0.75,0.84)
BraTS-T1CE	Brain Tumor	0.66 (0.53,0.80)	0.67 (0.50,0.80)	0.66 (0.54,0.79)	0.67 (0.52,0.80)	0.60 (0.49,0.76)	0.28 (0.20,0.40)	0.47 (0.34,0.62)	0.82 (0.60,0.90)	0.86 (0.66,0.91)	0.86 (0.66,0.91)	0.80 (0.55,0.89)
ISLES2022-DWI	Ischemic Stroke	0.56 (0.33,0.83)	0.58 (0.28,0.75)	0.63 (0.47,0.78)	0.54 (0.22,0.78)	0.60 (0.43,0.73)	0.06 (0.03,0.18)	0.73 (0.65,0.78)	0.81 (0.61,0.86)	0.77 (0.60,0.86)	0.80 (0.64,0.88)	0.80 (0.70,0.88)
ISLES2022-ADC	Ischemic Stroke	0.40 (0.14,0.59)	0.33 (0.13,0.60)	0.43 (0.18,0.59)	0.33 (0.12,0.61)	0.37 (0.30,0.52)	0.07 (0.03,0.15)	0.60 (0.56,0.67)	0.51 (0.42,0.71)	0.58 (0.38,0.72)	0.55 (0.42,0.73)	0.59 (0.45,0.71)
WMH-T1	White Matter Hyperintensity	0.07 (0.03,0.14)	0.08 (0.04,0.17)	0.13 (0.06,0.29)	0.09 (0.04,0.17)	0.38 (0.25,0.48)	0.02 (0.01,0.04)	0.38 (0.29,0.50)	0.48 (0.32,0.61)	0.43 (0.25,0.61)	0.53 (0.40,0.64)	0.54 (0.40,0.65)
WMH-FLAIR	White Matter Hyperintensity	0.13 (0.04,0.29)	0.11 (0.04,0.24)	0.17 (0.08,0.44)	0.15 (0.05,0.31)	0.42 (0.25,0.53)	0.02 (0.01,0.06)	0.45 (0.36,0.55)	0.59 (0.39,0.71)	0.46 (0.25,0.66)	0.58 (0.38,0.73)	0.54 (0.41,0.72)
PI-CAI	Prostate	0.81 (0.72,0.86)	0.78 (0.71,0.85)	0.78 (0.72,0.83)	0.77 (0.68,0.83)	0.70 (0.54,0.78)	0.24 (0.16,0.34)	0.88 (0.86,0.89)	0.92 (0.90,0.94)	0.91 (0.89,0.93)	0.92 (0.89,0.94)	0.91 (0.89,0.94)
Prostate-T2	Prostate	0.89 (0.86,0.92)	0.89 (0.85,0.91)	0.88 (0.85,0.90)	0.88 (0.85,0.91)	0.77 (0.60,0.83)	0.37 (0.28,0.49)	0.89 (0.88,0.90)	0.93 (0.91,0.94)	0.93 (0.91,0.94)	0.93 (0.90,0.94)	0.93 (0.91,0.94)
Prostate-ADC	Prostate	0.83 (0.81,0.88)	0.85 (0.82,0.88)	0.85 (0.81,0.88)	0.86 (0.82,0.88)	0.78 (0.66,0.84)	0.76 (0.72,0.80)	0.89 (0.87,0.90)	0.93 (0.91,0.94)	0.93 (0.91,0.94)	0.92 (0.90,0.94)	0.92 (0.90,0.93)
QIN-Prostate	Prostate	0.86 (0.82,0.89)	0.89 (0.86,0.90)	0.84 (0.80,0.90)	0.86 (0.82,0.90)	0.81 (0.73,0.86)	0.25 (0.23,0.53)	0.87 (0.79,0.89)	0.92 (0.90,0.93)	0.92 (0.91,0.93)	0.91 (0.90,0.92)	0.91 (0.91,0.92)
PROMISE12	Prostate	0.72 (0.66,0.76)	0.64 (0.57,0.71)	0.60 (0.48,0.70)	0.60 (0.44,0.65)	0.17 (0.11,0.21)	0.28 (0.19,0.37)	0.87 (0.84,0.89)	0.89 (0.87,0.90)	0.86 (0.84,0.88)	0.88 (0.87,0.89)	0.88 (0.87,0.90)
MSD-Prostate	Prostate Central Zone	0.58 (0.57,0.67)	0.58 (0.54,0.64)	0.60 (0.59,0.63)	0.59 (0.53,0.64)	0.27 (0.19,0.39)	0.54 (0.51,0.58)	0.87 (0.86,0.88)	0.75 (0.74,0.77)	0.74 (0.71,0.78)	0.69 (0.68,0.76)	0.73 (0.72,0.75)
MSD-Prostate	Prostate Peripheral Zone	0.19 (0.15,0.45)	0.20 (0.14,0.46)	0.26 (0.15,0.46)	0.24 (0.11,0.45)	0.26 (0.24,0.42)	0.25 (0.22,0.32)	0.51 (0.49,0.61)	0.35 (0.23,0.52)	0.28 (0.25,0.60)	0.29 (0.18,0.55)	0.35 (0.29,0.62)
QIN-Lesion	Prostate Peripheral Zone	0.65 (0.55,0.70)	0.62 (0.46,0.73)	0.54 (0.51,0.65)	0.62 (0.52,0.65)	0.60 (0.59,0.73)	0.35 (0.22,0.43)	0.62 (0.40,0.71)	0.70 (0.53,0.78)	0.68 (0.54,0.69)	0.65 (0.53,0.77)	0.66 (0.55,0.74)
ZIB-OAI	Bone(Femur)	0.95 (0.90,0.97)	0.96 (0.94,0.97)	0.97 (0.94,0.97)	0.96 (0.92,0.97)	0.45 (0.02,0.71)	0.77 (0.72,0.80)	0.75 (0.72,0.78)	0.98 (0.97,0.98)	0.97 (0.97,0.98)	0.97 (0.97,0.98)	0.98 (0.97,0.98)
ZIB-OAI	Bone(Tibia)	0.97 (0.96,0.98)	0.95 (0.94,0.96)	0.96 (0.96,0.97)	0.96 (0.93,0.97)	0.33 (0.02,0.54)	0.78 (0.74,0.82)	0.82 (0.80,0.83)	0.98 (0.98,0.98)	0.96 (0.95,0.98)	0.98 (0.97,0.98)	0.98 (0.97,0.98)
ZIB-OAI	Cartilage(Femur)	0.64 (0.55,0.69)	0.63 (0.57,0.66)	0.57 (0.47,0.66)	0.61 (0.51,0.66)	0.02 (0.00,0.05)	0.00 (0.00,0.00)	0.21 (0.19,0.23)	0.85 (0.80,0.86)	0.85 (0.83,0.87)	0.85 (0.82,0.87)	0.83 (0.81,0.86)
ZIB-OAI	Cartilage(Tibia)	0.47 (0.39,0.56)	0.42 (0.36,0.46)	0.39 (0.31,0.46)	0.37 (0.31,0.46)	0.22 (0.00,0.41)	0.19 (0.05,0.32)	0.51 (0.47,0.54)	0.72 (0.64,0.78)	0.72 (0.67,0.77)	0.74 (0.71,0.77)	0.70 (0.65,0.77)
OAI-iMorphics-DESS	Bone(Femur)	0.93 (0.87,0.95)	0.95 (0.93,0.96)	0.95 (0.58,0.97)	0.94 (0.90,0.96)	0.62 (0.02,0.81)	0.56 (0.43,0.74)	0.77 (0.75,0.80)	0.97 (0.97,0.98)	0.98 (0.97,0.98)	0.97 (0.97,0.98)	0.97 (0.97,0.98)
OAI-iMorphics-DESS	Cartilage(Femur)	0.61 (0.56,0.68)	0.56 (0.50,0.59)	0.44 (0.12,0.57)	0.53 (0.43,0.61)	0.00 (0.00,0.01)	0.00 (0.00,0.03)	0.12 (0.11,0.15)	0.82 (0.77,0.84)	0.82 (0.80,0.85)	0.81 (0.80,0.84)	0.82 (0.78,0.84)
OAI-iMorphics-DESS	Bone(Tibia)	0.97 (0.96,0.97)	0.96 (0.95,0.97)	0.97 (0.96,0.97)	0.96 (0.95,0.97)	0.13 (0.02,0.52)	0.43 (0.28,0.53)	0.82 (0.77,0.84)	0.98 (0.98,0.98)	0.97 (0.95,0.98)	0.98 (0.98,0.98)	0.98 (0.98,0.98)
OAI-iMorphics-DESS	Cartilage(Tibia)	0.48 (0.31,0.54)	0.39 (0.31,0.46)	0.34 (0.27,0.42)	0.37 (0.29,0.45)	0.25 (0.01,0.42)	0.23 (0.13,0.32)	0.52 (0.48,0.54)	0.76 (0.68,0.80)	0.76 (0.70,0.79)	0.77 (0.73,0.79)	0.74 (0.68,0.78)
OAI-iMorphics-DESS	Patella	0.63 (0.53,0.78)	0.92 (0.86,0.95)	0.96 (0.95,0.96)	0.89 (0.60,0.95)	0.33 (0.04,0.56)	0.70 (0.55,0.76)	0.82 (0.82,0.84)	0.96 (0.95,0.97)	0.96 (0.94,0.97)	0.96 (0.95,0.97)	0.95 (0.92,0.96)
OAI-iMorphics-DESS	Cartilage(Patella)	0.32 (0.25,0.48)	0.49 (0.29,0.63)	0.40 (0.23,0.66)	0.29 (0.24,0.61)	0.04 (0.03,0.15)	0.24 (0.17,0.31)	0.53 (0.47,0.60)	0.84 (0.74,0.87)	0.84 (0.74,0.87)	0.83 (0.77,0.86)	0.83 (0.77,0.86)
OAI-iMorphics-TSE	Bone(Femur)	0.58 (0.54,0.62)	0.87 (0.61,0.92)	0.50 (0.36,0.55)	0.57 (0.52,0.65)	0.87 (0.81,0.90)	0.73 (0.60,0.81)	0.81 (0.78,0.82)	0.93 (0.64,0.94)	0.93 (0.63,0.94)	0.92 (0.74,0.95)	0.86 (0.69,0.93)
OAI-iMorphics-TSE	Cartilage(Femur)	0.07 (0.05,0.14)	0.47 (0.12,0.55)	0.04 (0.02,0.06)	0.07 (0.05,0.16)	0.03 (0.03,0.04)	0.05 (0.03,0.08)	0.10 (0.08,0.11)	0.44 (0.11,0.62)	0.55 (0.15,0.67)	0.53 (0.14,0.67)	0.32 (0.13,0.62)
OAI-iMorphics-TSE	Bone(Tibia)	0.92 (0.91,0.93)	0.91 (0.89,0.93)	0.92 (0.90,0.93)	0.84 (0.79,0.88)	0.73 (0.64,0.84)	0.78 (0.69,0.84)	0.83 (0.76,0.85)	0.92 (0.91,0.94)	0.91 (0.89,0.93)	0.92 (0.91,0.94)	0.92 (0.90,0.94)
OAI-iMorphics-TSE	Cartilage(Tibia)	0.39 (0.32,0.44)	0.36 (0.33,0.41)	0.30 (0.26,0.33)	0.31 (0.27,0.34)	0.34 (0.30,0.41)	0.15 (0.10,0.19)	0.47 (0.43,0.48)	0.56 (0.53,0.59)	0.56 (0.53,0.59)	0.57 (0.46,0.60)	0.53 (0.46,0.59)
OAI-iMorphics-TSE	Patella	0.77 (0.61,0.86)	0.84 (0.75,0.89)	0.88 (0.80,0.90)	0.60 (0.42,0.78)	0.71 (0.65,0.80)	0.60 (0.54,0.68)	0.81 (0.77,0.82)	0.89 (0.86,0.91)	0.89 (0.84,0.91)	0.90 (0.86,0.92)	0.91 (0.84,0.92)
OAI-iMorphics-TSE	Cartilage(Patella)	0.39 (0.30,0.47)	0.36 (0.28,0.43)	0.26 (0.19,0.43)	0.23 (0.05,0.38)	0.23 (0.05,0.38)	0.48 (0.42,0.52)	0.67 (0.57,0.72)	0.86 (0.75,0.72)	0.86 (0.75,0.72)	0.64 (0.54,0.70)	0.63 (0.56,0.71)
OAI-iMorphics-FLASH	Bone(Femur)	0.93 (0.86,0.94)	0.92 (0.91,0.95)	0.96 (0.80,0.97)	0.94 (0.90,0.96)	0.94 (0.93,0.94)	0.43 (0.38,0.49)	0.73 (0.73,0.75)	0.97 (0.97,0.97)	0.97 (0.96,0.97)	0.97 (0.95,0.97)	0.97 (0.96,0.97)
OAI-iMorphics-FLASH	Cartilage(Femur)	0.53 (0.51,0.56)	0.53 (0.49,0.56)	0.57 (0.25,0.62)	0.51 (0.44,0.55)	0.01 (0.00,0.01)	0.04 (0.02,0.14)	0.15 (0.13,0.17)	0.73 (0.70,0.77)	0.76 (0.75,0.77)	0.76 (0.75,0.77)	0.74 (0.71,0.75)
OAI-iMorphics-FLASH	Bone(Tibia)	0.92 (0.55,0.95)	0.94 (0.93,0.95)	0.96 (0.95,0.96)	0.95 (0.92,0.96)	0.88 (0.80,0.92)	0.42 (0.14,0.64)	0.85 (0.84,0.87)	0.97 (0.97,0.98)	0.92 (0.63,0.93)	0.97 (0.94,0.97)	0.97 (0.97,0.98)
OAI-iMorphics-FLASH	Cartilage(Tibia)	0.50 (0.39,0.59)	0.48 (0.44,0.53)	0.48 (0.42,0.54)	0.38 (0.29,0.40)	0.55 (0.50,0.59)	0.22 (0.11,0.35)	0.50 (0.48,0.51)	0.71 (0.65,0.73)	0.69 (0.63,0.72)	0.71 (0.68,0.73)	0.69 (0.65,0.73)
OAI-iMorphics-FLASH	Patella	0.87 (0.79,0.90)	0.90 (0.74,0.93)	0.93 (0.93,0.94)	0.87 (0.83,0.89)	0.85 (0.81,0.86)	0.76 (0.68,0.80)	0.82 (0.80,0.84)	0.94 (0.93,0.95)	0.93 (0.89,0.94)	0.93 (0.92,0.94)	0.93 (0.90,0.94)
OAI-iMorphics-FLASH	Cartilage(Patella)	0.31 (0.19,0.45)	0.34 (0.27,0.45)	0.29 (0.27,0.53)	0.34 (0.19,0.42)	0.55 (0.52,0.61)	0.25 (0.12,0.38)	0.56 (0.54,0.61)	0.64 (0.56,0.74)	0.61 (0.61,0.76)	0.64 (0.62,0.76)	0.68 (0.62,0.76)
OAI-AKOA	Bone(Femur)	0.89 (0.88,0.90)	0.91 (0.91,0.92)	0.70 (0.60,0.80)	0.68 (0.58,0.77)	0.62 (0.80,0.84)	0.76 (0.75,0.78)	0.79 (0.79,0.79)	0.96 (0.95,0.97)	0.96 (0.96,0.97)	0.96 (0.95,0.97)	0.96 (0.94,0.97)
OAI-AKOA	Bone(Tibia)	0.97 (0.97,0.97)	0.95 (0.94,0.95)	0.97 (0.96,0.97)	0.95 (0.94,0.96)	0.85 (0.81,0.82)	0.82 (0.81,0.82)	0.84 (0.84,0.84)	0.98 (0.98,0.98)	0.96 (0.96,0.97)	0.98 (0.98,0.98)	0.98 (0.98,0.98)
OAI-AKOA	Cartilage(Femur)	0.64 (0.61,0.66)	0.47 (0.43,0.51)	0.29 (0.15,0.44)	0.13 (0.07,0.20)	0.01 (0.01,0.01)	0.00 (0.00,0.00)	0.13 (0.13,0.14)	0.77 (0.74,0.80)	0.77 (0.74,0.80)	0.77 (0.74,0.79)	0.75 (0.72,0.78)
OAI-AKOA	Cartilage(Patella)	0.40 (0.38,0.41)	0.37 (0.36,0.37)	0.35 (0.29,0.40)	0.53 (0.52,0.54)	0.14 (0.08,0.20)	0.13 (0.13,0.13)	0.63 (0.61,0.65)	0.73 (0.72,0.74)	0.74 (0.71,0.77)	0.75 (0.73,0.76)	0.71 (0.67,0.75)
OAI-AKOA	Cartilage(Tibia)	0.48 (0.44,0.53)	0.37 (0.37,0.38)	0.43 (0.41,0.45)	0.59 (0.56,0.63)	0.46 (0.42,0.50)	0.12 (0.08,0.17)	0.54 (0.52,0.56)	0.73 (0.69,0.77)	0.75 (0.71,0.78)	0.74 (0.71,0.77)	0.72 (0.69,0.75)
OAI-AKOA	Patella	0.55 (0.52,0.59)	0.86 (0.86,0.87)	0.88 (0.84,0.91)	0.75 (0.72,0.79)	0.40 (0.22,0.57)	0.86 (0.86,0.87)	0.82 (0.81,0.83)	0.95 (0.95,0.95)	0.95 (0.95,0.95)	0.94 (0.94,0.94)	0.94 (0.94,0.94)
MSK-T2	Bone(Femur)	0.94 (0.94,0.95)	0.94 (0.94,0.94)	0.95 (0.94,0.95)	0.94 (0.93,0.95)	0.75 (0.73,0.77)	0.67 (0.61,0.72)	0.76 (0.70,0.81)	0.96 (0.96,0.96)	0.96 (0.96,0.96)	0.96 (0.95,0.96)	0.96 (0.95,0.96)
MSK-T2	Bone(Tibia)	0.95 (0.92,0.96)	0.95 (0.94,0.96)	0.95 (0.90,0.96)	0.95 (0.90,0.96)	0.26 (0.01,0.57)	0.87 (0.85,0.88)	0.86 (0.84,0.86)	0.96 (0.95,0.96)	0.93 (0.92,0.94)	0.93 (0.91,0.95)	0.95 (0.94,0.95)
MSK-T2	Cartilage(Femur)	0.70 (0.67,0.73)	0.66 (0.66,0.66)	0.77 (0.75,0.78)	0.75 (0.75,0.75)	0.03 (0.02,0.04)	0.01 (0.01,0.07)	0.08 (0.07,0.09)	0.87 (0.85,0.88)	0.87 (0.85,0.88)	0.87 (0.86,0.88)	0.87 (0.85,0.88)
MSK-T2	Cartilage(Patella)	0.84 (0.83,0.87)	0.82 (0.79,0.85)	0.83 (0.79,0.86)	0.70 (0.53,0.74)	0.48 (0.40,0.51)	0.15 (0.10,0.22)	0.70 (0.69,0.72)	0.89 (0.88,0.90)	0.89 (0.88,0.90)	0.89 (0.88,0.90)	0.89 (0.88,0.90)
MSK-T2	Cartilage(Tibia)	0.40 (0.36,0.46)	0.31 (0.28,0.37)	0.45 (0.37,0.56)	0.39 (0.36,0.41)	0.33 (0.20,0.43)	0.33 (0.21,0.39)	0.58 (0.57,0.60)	0.61 (0.57,0.63)	0.38 (0.32,0.49)	0.68 (0.55,0.77)	0.55 (0.49,0.60)
MSK-T2	Patella	0.94 (0.92,0.94)	0.92 (0.89,0.93)	0.93 (0.91,0.94)	0.90 (0.76,0.92)	0.78 (0.67,0.86)	0.83 (0.82,0.85)	0.93 (0.93,0.94)	0.95 (0.93,0.95)	0.94 (0.92,0.95)	0.94 (0.92,0.95)	0.94 (0.92,0.94)
MSK-FLASH	Bone(Femur)	0.03 (0.03,0.04)</										

Dataset	Target	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2dloss	GVT-v1	GVT-v2	SAMRI-3D
AMOS-MR	Aorta	0.87 (0.59,0.90)	0.87 (0.35,0.89)	0.87 (0.60,0.89)	0.83 (0.28,0.89)	0.69 (0.48,0.78)	0.41 (0.39,0.48)	0.80 (0.74,0.84)	0.85 (0.63,0.90)	0.89 (0.77,0.91)	0.89 (0.85,0.90)	0.89 (0.31,0.90)
AMOS-MR	Duodenum	0.13 (0.07,0.19)	0.14 (0.12,0.28)	0.23 (0.22,0.27)	0.18 (0.16,0.24)	0.26 (0.16,0.30)	0.14 (0.05,0.26)	0.54 (0.41,0.60)	0.37 (0.31,0.44)	0.33 (0.24,0.48)	0.31 (0.18,0.46)	0.42 (0.35,0.50)
AMOS-MR	Esophagus	0.43 (0.11,0.50)	0.23 (0.10,0.37)	0.24 (0.05,0.36)	0.08 (0.05,0.18)	0.34 (0.24,0.40)	0.26 (0.24,0.28)	0.62 (0.61,0.64)	0.57 (0.42,0.69)	0.43 (0.38,0.49)	0.55 (0.46,0.70)	0.57 (0.49,0.66)
AMOS-MR	Gallbladder	0.16 (0.03,0.40)	0.05 (0.01,0.22)	0.17 (0.07,0.27)	0.11 (0.03,0.43)	0.53 (0.39,0.65)	0.43 (0.28,0.54)	0.76 (0.72,0.80)	0.31 (0.13,0.48)	0.24 (0.06,0.40)	0.46 (0.24,0.55)	0.45 (0.20,0.64)
AMOS-MR	Inferior Vena Cava	0.31 (0.22,0.35)	0.44 (0.20,0.67)	0.32 (0.20,0.58)	0.36 (0.20,0.57)	0.28 (0.25,0.37)	0.36 (0.28,0.42)	0.76 (0.70,0.78)	0.69 (0.49,0.77)	0.67 (0.55,0.79)	0.65 (0.57,0.74)	0.70 (0.65,0.77)
AMOS-MR	Left Adrenal Gland	0.14 (0.05,0.32)	0.10 (0.07,0.12)	0.08 (0.05,0.14)	0.10 (0.06,0.12)	0.20 (0.06,0.25)	0.15 (0.13,0.19)	0.47 (0.43,0.51)	0.49 (0.31,0.66)	0.52 (0.37,0.64)	0.52 (0.40,0.63)	0.53 (0.35,0.66)
AMOS-MR	Left Kidney	0.68 (0.63,0.78)	0.64 (0.51,0.67)	0.81 (0.80,0.84)	0.74 (0.71,0.80)	0.73 (0.62,0.79)	0.79 (0.72,0.83)	0.85 (0.82,0.86)	0.91 (0.89,0.93)	0.90 (0.88,0.92)	0.91 (0.83,0.92)	0.90 (0.86,0.91)
AMOS-MR	Liver	0.90 (0.83,0.94)	0.87 (0.86,0.88)	0.90 (0.87,0.93)	0.87 (0.75,0.91)	0.65 (0.58,0.73)	0.17 (0.04,0.30)	0.79 (0.77,0.82)	0.91 (0.86,0.94)	0.90 (0.84,0.94)	0.88 (0.84,0.93)	0.92 (0.86,0.94)
AMOS-MR	Pancreas	0.28 (0.16,0.42)	0.28 (0.18,0.40)	0.29 (0.21,0.46)	0.30 (0.07,0.63)	0.32 (0.24,0.36)	0.11 (0.05,0.26)	0.67 (0.60,0.70)	0.72 (0.61,0.75)	0.58 (0.29,0.73)	0.72 (0.57,0.74)	0.74 (0.61,0.75)
AMOS-MR	Right Adrenal Gland	0.02 (0.01,0.05)	0.03 (0.02,0.19)	0.02 (0.02,0.04)	0.02 (0.02,0.05)	0.31 (0.26,0.41)	0.10 (0.06,0.14)	0.43 (0.39,0.49)	0.40 (0.31,0.42)	0.36 (0.34,0.38)	0.46 (0.36,0.51)	0.44 (0.37,0.48)
AMOS-MR	Right Kidney	0.00 (0.00,0.36)	0.00 (0.00,0.00)	0.00 (0.00,0.01)	0.00 (0.00,0.02)	0.73 (0.46,0.78)	0.78 (0.58,0.81)	0.82 (0.78,0.85)	0.90 (0.89,0.90)	0.32 (0.01,0.82)	0.90 (0.89,0.91)	0.90 (0.83,0.91)
AMOS-MR	Spleen	0.56 (0.47,0.64)	0.52 (0.44,0.56)	0.83 (0.78,0.88)	0.79 (0.76,0.86)	0.61 (0.52,0.68)	0.61 (0.41,0.77)	0.83 (0.78,0.86)	0.93 (0.89,0.94)	0.92 (0.88,0.93)	0.87 (0.78,0.91)	0.90 (0.84,0.94)
AMOS-MR	Stomach	0.67 (0.48,0.80)	0.51 (0.42,0.60)	0.50 (0.45,0.67)	0.59 (0.40,0.71)	0.42 (0.41,0.43)	0.06 (0.02,0.08)	0.74 (0.72,0.77)	0.77 (0.68,0.79)	0.73 (0.67,0.83)	0.77 (0.60,0.83)	0.65 (0.50,0.79)
MSK-Shoulder	Humerus cartilage	0.45 (0.36,0.56)	0.43 (0.29,0.49)	0.43 (0.33,0.57)	0.39 (0.32,0.40)	0.39 (0.39,0.44)	0.30 (0.22,0.36)	0.58 (0.49,0.59)	0.74 (0.67,0.76)	0.64 (0.62,0.71)	0.76 (0.73,0.78)	0.76 (0.73,0.79)
MSK-Shoulder	Scapula cartilage	0.00 (0.00,0.00)	0.01 (0.00,0.47)	0.00 (0.00,0.00)	0.00 (0.00,0.62)	0.02 (0.02,0.04)	0.03 (0.01,0.07)	0.14 (0.11,0.15)	0.05 (0.05,0.31)	0.56 (0.04,0.77)	0.05 (0.01,0.07)	0.05 (0.03,0.14)
CervicalCancer	Cervix Cancer	0.70 (0.70,0.70)	0.67 (0.67,0.67)	0.63 (0.63,0.63)	0.66 (0.66,0.66)	0.53 (0.53,0.53)	0.55 (0.55,0.55)	0.68 (0.68,0.68)	0.74 (0.74,0.74)	0.75 (0.75,0.75)	0.73 (0.73,0.73)	0.76 (0.76,0.76)
Mean $\pm$ SD		0.57 $\pm$ 0.29	0.56 $\pm$ 0.30	0.56 $\pm$ 0.30	0.54 $\pm$ 0.30	0.45 $\pm$ 0.26	0.39 $\pm$ 0.28	0.63 $\pm$ 0.22	0.74 $\pm$ 0.23	0.73 $\pm$ 0.23	0.77 $\pm$ 0.20	0.76 $\pm$ 0.20

Table S10. HD by target — median (Q1, Q3) across volumes; final row is mean $\pm$ SD.

Dataset	Target	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2loss	GVT-v1	GVT-v2	SAMRI-3D
MSD-BrainTumour	Edema	58.89 (49.18,68.89)	54.80 (43.02,66.20)	54.90 (43.38,65.55)	61.06 (48.31,72.89)	48.49 (36.74,57.79)	23.26 (15.66,31.08)	19.25 (14.27,24.19)	42.28 (29.72,52.45)	32.65 (21.37,46.05)	49.13 (31.14,74.48)	42.20 (27.20,52.34)
MSD-BrainTumour	Enhancing tumour	68.17 (57.78,78.93)	63.11 (45.55,75.98)	61.26 (37.55,73.72)	70.81 (56.02,80.91)	45.77 (26.41,65.17)	20.86 (13.83,28.28)	19.30 (7.86,17.27)	27.25 (18.05,45.15)	19.30 (11.67,30.59)	31.96 (18.99,58.14)	27.20 (16.05,53.46)
MSD-BrainTumour	Non-enhancing tumour	65.68 (55.75,76.42)	61.12 (41.43,72.63)	57.12 (34.41,74.30)	65.06 (45.02,78.22)	49.49 (32.12,65.43)	21.15 (13.99,29.38)	11.87 (8.06,20.62)	34.43 (22.05,49.00)	30.81 (21.25,45.07)	39.50 (26.01,64.78)	38.48 (23.63,53.80)
BraTS-FLAIR	Brain Tumor	22.19 (14.04,30.95)	21.42 (15.69,37.59)	21.54 (11.79,34.70)	24.93 (16.06,40.72)	13.96 (10.92,19.40)	13.96 (10.92,19.40)	19.44 (11.57,53.37)	14.05 (8.77,24.62)	14.92 (9.68,22.60)	13.25 (8.46,25.81)	13.38 (8.53,21.85)
BraTS-T1	Brain Tumor	24.71 (17.51,35.07)	29.91 (19.80,37.10)	29.05 (20.47,38.46)	34.43 (26.21,51.98)	27.08 (23.33,34.02)	20.00 (15.57,26.21)	17.01 (11.97,21.45)	18.07 (13.42,26.81)	18.03 (14.40,27.24)	19.40 (12.61,31.70)	16.76 (12.78,23.19)
BraTS-T1CE	Brain Tumor	8.80 (6.32,11.58)	9.27 (6.65,13.38)	8.80 (6.71,12.44)	9.77 (7.35,16.26)	11.18 (9.43,15.20)	11.18 (9.43,15.20)	8.15 (6.16,12.53)	7.62 (5.51,10.68)	7.07 (5.20,10.63)	7.45 (5.02,10.64)	7.75 (5.54,11.52)
ISLES2022-DWI	Ischemic Stroke	16.54 (8.47,26.74)	17.42 (10.46,25.51)	11.77 (4.54,26.05)	19.10 (9.50,26.94)	5.61 (2.24,9.84)	16.03 (13.03,19.12)	2.83 (2.00,4.21)	3.74 (1.49,7.45)	3.87 (2.06,7.88)	3.24 (1.41,5.17)	3.46 (1.80,5.98)
ISLES2022-ADC	Ischemic Stroke	16.58 (5.10,26.60)	20.55 (8.00,24.57)	10.03 (5.46,18.69)	22.13 (8.97,29.00)	10.03 (5.49,17.22)	16.36 (14.01,19.67)	3.16 (2.40,5.02)	4.69 (2.93,8.65)	4.58 (3.12,8.80)	4.47 (2.96,7.98)	4.24 (2.45,7.87)
WMH-T1	White Matter Hyperintensity	45.84 (32.75,62.95)	44.37 (32.83,66.59)	33.84 (14.99,54.95)	50.98 (33.83,68.88)	7.00 (3.32,16.90)	16.96 (14.12,22.62)	5.15 (3.61,11.12)	6.78 (3.46,11.26)	8.60 (3.94,14.11)	4.53 (2.40,8.08)	4.58 (2.45,9.26)
WMH-FLAIR	White Matter Hyperintensity	55.18 (33.77,83.63)	54.84 (35.56,82.97)	48.21 (23.17,74.44)	56.87 (37.43,84.24)	8.18 (3.74,20.48)	16.90 (14.20,21.12)	6.89 (3.61,53.31)	7.78 (4.12,17.67)	15.23 (5.81,27.65)	6.56 (3.04,16.75)	6.00 (3.32,16.20)
PI-CAI	Prostate	34.07 (15.61,53.84)	40.14 (24.41,60.61)	41.78 (22.73,69.01)	56.59 (34.00,87.18)	45.11 (30.19,60.59)	40.52 (34.21,50.68)	14.63 (11.67,17.66)	6.56 (4.69,8.86)	8.15 (5.83,11.69)	8.03 (5.05,13.81)	8.22 (5.17,13.18)
Prostate-T2	Prostate	11.38 (7.26,15.51)	11.46 (7.03,18.71)	13.15 (8.11,17.90)	14.88 (8.87,21.41)	23.16 (13.71,40.95)	35.90 (25.53,44.73)	13.15 (10.51,15.30)	6.00 (5.02,7.12)	6.04 (5.02,7.14)	6.12 (5.12,8.02)	5.92 (5.00,7.07)
Prostate-ADC	Prostate	9.38 (5.10,12.33)	6.16 (5.10,12.33)	8.60 (4.58,10.82)	7.28 (4.47,12.08)	8.72 (5.92,12.53)	7.21 (6.16,8.60)	2.83 (2.45,3.61)	2.45 (2.24,3.00)	2.83 (2.24,3.46)	2.83 (2.24,3.46)	2.83 (2.24,3.61)
QIN-Prostate	Prostate	18.06 (10.00,32.16)	11.18 (9.49,18.79)	24.54 (10.25,42.56)	13.34 (11.45,33.66)	24.56 (18.14,29.53)	37.54 (22.16,70.72)	11.58 (10.63,14.32)	7.14 (7.00,9.90)	9.49 (6.40,10.49)	8.60 (7.68,20.10)	8.79 (6.48,9.95)
PROMISE12	Prostate	47.98 (42.36,87.95)	50.79 (48.72,91.99)	86.98 (56.81,116.98)	86.27 (69.00,100.63)	35.83 (11.54,43.55)	25.83 (22.24,54.64)	10.77 (6.52,14.23)	6.48 (4.74,9.72)	10.00 (8.71,17.49)	5.92 (4.79,8.20)	7.07 (6.08,9.35)
MSD-Prostate	Prostate Central Zone	52.07 (48.92,57.97)	59.04 (33.90,56.14)	45.57 (33.90,56.14)	40.39 (28.38,47.92)	29.04 (24.86,51.02)	40.39 (24.86,51.02)	17.87 (14.99,20.23)	7.87 (7.18,8.29)	9.15 (7.59,10.47)	8.68 (7.93,9.69)	7.98 (7.35,9.65)
MSD-Prostate	Prostate Peripheral Zone	38.95 (18.70,55.42)	39.33 (16.75,49.53)	42.16 (18.22,58.95)	39.10 (15.82,42.54)	42.59 (36.25,53.76)	21.55 (18.37,23.67)	13.04 (10.89,13.93)	14.61 (11.35,18.53)	12.66 (10.86,14.23)	15.90 (8.92,20.96)	12.77 (10.38,17.85)
QIN-Lesion	Prostate Peripheral Zone	9.43 (5.73,17.18)	4.00 (4.03,12.04)	8.50 (4.53,12.26)	7.35 (4.87,10.46)	8.60 (4.96,11.52)	12.57 (11.52,13.74)	7.87 (5.34,10.00)	5.92 (3.86,10.23)	5.22 (3.99,10.21)	5.74 (4.30,9.51)	7.87 (2.99,11.25)
ZIB-OAI	Bone(Femur)	39.52 (20.64,91.34)	37.72 (18.55,64.29)	27.16 (14.26,51.31)	60.04 (25.53,117.73)	71.63 (61.06,90.45)	39.87 (34.14,45.59)	54.27 (50.45,58.61)	13.23 (9.81,20.41)	12.59 (9.21,18.32)	11.70 (8.50,18.62)	17.86 (11.83,41.69)
ZIB-OAI	Bone(Tibia)	52.93 (11.30,113.39)	68.10 (59.66,85.53)	66.19 (56.26,71.82)	90.00 (63.83,126.85)	60.19 (57.99,80.67)	40.19 (35.04,43.95)	29.03 (26.87,32.74)	9.11 (7.11,12.74)	55.94 (10.16,10.19)	10.22 (7.25,10.22)	20.04 (10.24,53.46)
ZIB-OAI	Cartilage(Femur)	137.31 (90.60,155.84)	120.66 (84.62,146.99)	114.04 (82.05,139.44)	126.42 (97.32,154.23)	90.19 (84.98,103.14)	86.66 (83.75,88.80)	75.34 (70.70,79.85)	17.80 (13.71,28.54)	19.26 (14.04,27.09)	16.40 (13.37,21.00)	19.85 (15.07,24.49)
ZIB-OAI	Cartilage(Tibia)	85.42 (70.07,121.55)	117.61 (93.17,174.71)	109.58 (93.11,154.22)	116.56 (101.02,153.62)	73.74 (55.88,87.12)	51.38 (42.86,57.82)	47.57 (39.22,58.81)	19.76 (15.82,26.77)	18.75 (16.82,26.77)	18.87 (15.84,25.26)	21.57 (18.31,32.08)
OAI-iMorphics-DESS	Bone(Femur)	43.02 (24.43,97.00)	35.15 (23.90,55.78)	41.14 (26.54,68.86)	67.98 (37.00,119.32)	80.50 (59.47,91.62)	53.11 (47.79,62.95)	54.36 (51.44,57.09)	14.11 (9.79,21.12)	15.20 (9.97,20.11)	13.86 (10.23,17.84)	20.35 (10.47,31.84)
OAI-iMorphics-DESS	Cartilage(Femur)	125.32 (81.74,188.76)	115.44 (98.78,141.37)	98.25 (87.23,127.12)	149.67 (96.05,171.62)	87.60 (84.70,106.43)	78.80 (73.85,83.98)	74.76 (72.24,79.94)	19.62 (15.65,32.52)	18.35 (15.54,27.26)	16.81 (14.84,29.79)	21.42 (16.58,26.65)
OAI-iMorphics-DESS	Bone(Tibia)	35.75 (11.42,99.92)	61.11 (37.52,71.00)	65.51 (54.76,72.82)	86.17 (71.20,118.13)	62.20 (56.58,68.50)	67.02 (55.28,74.12)	30.25 (27.53,33.43)	7.55 (6.40,11.05)	16.75 (9.31,62.04)	9.51 (7.05,12.54)	10.31 (7.00,29.54)
OAI-iMorphics-DESS	Cartilage(Tibia)	110.74 (78.51,149.99)	150.24 (104.28,187.92)	132.24 (106.12,166.93)	114.73 (94.18,138.35)	76.38 (53.62,82.15)	39.27 (32.20,43.57)	52.15 (31.28,57.58)	21.43 (16.96,26.63)	21.52 (19.43,25.46)	20.01 (17.20,24.96)	21.34 (19.23,26.55)
OAI-iMorphics-DESS	Patella	225.97 (198.96,256.61)	161.14 (170.228,55)	9.70 (7.05,14.96)	172.90 (77.45,275.62)	61.43 (45.96,82.09)	36.06 (33.34,38.67)	15.37 (13.40,18.74)	8.94 (7.23,11.44)	10.05 (7.00,74.23)	9.08 (7.35,27.37)	29.94 (27.90,40.40)
OAI-iMorphics-DESS	Cartilage(Patella)	196.71 (166.03,226.07)	168.58 (139.35,244.80)	146.47 (98.83,169.41)	171.93 (137.58,199.13)	51.28 (46.47,56.75)	60.73 (57.80,64.01)	12.31 (10.76,15.28)	9.53 (6.38,13.82)	9.95 (6.44,14.07)	8.97 (5.72,12.67)	8.79 (6.04,12.29)
OAI-iMorphics-TSE	Bone(Femur)	46.59 (41.63,57.83)	45.52 (30.69,70.94)	49.52 (43.45,64.72)	53.21 (40.51,72.10)	37.84 (31.29,56.43)	51.32 (44.96,55.57)	50.31 (45.84,54.96)	21.06 (14.82,34.63)	22.02 (16.08,32.23)	25.73 (18.86,31.53)	27.34 (18.55,50.91)
OAI-iMorphics-TSE	Cartilage(Femur)	82.88 (71.63,97.49)	78.99 (67.91,133.58)	106.07 (79.39,122.69)	106.30 (81.96,119.97)	65.97 (57.63,86.53)	65.97 (57.63,86.53)	65.97 (64.81,71.86)	62.59 (56.38,77.08)	62.59 (56.38,77.08)	51.11 (38.61,67.21)	64.62 (40.18,73.40)
OAI-iMorphics-TSE	Bone(Tibia)	32.04 (22.21,101.33)	59.59 (70.14,112.79)	62.57 (41.83,81.26)	200.15 (136.12,223.99)	51.00 (43.02,55.57)	47.70 (39.14,50.83)	29.51 (26.43,34.69)	20.16 (17.82,18.78)	28.78 (24.72,37.59)	18.31 (14.76,23.29)	23.32 (18.41,29.03)
OAI-iMorphics-TSE	Cartilage(Tibia)	47.59 (36.88,82.56)	79.04 (36.44,99.78)	70.02 (52.42,88.60)	96.41 (83.23,132.32)	66.33 (40.74,78.62)	37.42 (32.26,45.54)	46.40 (40.25,45.54)	24.23 (17.61,44.36)	23.14 (17.38,35.98)	21.26 (15.69,40.83)	23.38 (18.37,34.61)
OAI-iMorphics-TSE	Patella	147.33 (95.55,160.37)	37.05 (23.79,83.07)	24.33 (17.85,59.50)	166.19 (41.53,216.59)	41.46 (27.49,49.83)	32.42 (26.18,41.74)	14.21 (11.91,22.23)	15.20 (11.42,22.06)	15.57 (11.48,24.85)	15.80 (10.76,21.84)	14.58 (10.58,21.52)
OAI-iMorphics-TSE	Cartilage(Patella)	90.10 (75.73,120.16)	86.39 (64.64,114.43)	100.56 (83.19,122.36)	150.90 (108.26,181.94)	41.65 (35.97,56.38)	56.42 (47.27,62.26)	13.47 (11.18,25.21)	17.36 (12.03,24.53)	16.42 (13.11,43.68)	16.21 (12.98,28.87)	15.58 (12.98,22.95)
OAI-iMorphics-FLASH	Bone(Femur)	107.13 (89.74,142.05)	73.23 (37.02,97.84)	36.57 (26.58,53.01)	101.09 (88.08,128.25)	33.38 (29.81,67.25)	68.77 (62.95,80.53)	66.22 (61.70,73.31)	12.77 (11.44,25.27)	29.77 (18.29,58.02)	37.11 (23.66,53.45)	37.11 (23.66,53.45)
OAI-iMorphics-FLASH	Cartilage(Femur)	162.34 (154.39,170.58)	153.94 (117.62,173.65)	89.81 (83.35,129.15)	145.66 (71.91,163.24)	92.28 (89.16,93.91)	73.71 (71.16,75.63)	35.74 (26.17,55.45)	26.51 (20.25,35.65)	26.64 (24.53,29.60)	30.06 (24.91,49.71)	30.06 (24.91,49.71)
OAI-iMorphics-FLASH	Bone(Tibia)	208.13 (143.33,261.16)	63.93 (52.39,87.77)	129.26 (35.12,234.30)	217.81 (136.57,250.34)	60.47 (54.61,66.08)	45.33 (33.05,54.01)	46.57 (35.00,55.79)	12.02 (9.85,14.09)	68.68 (62.53,229.92)	35.03 (10.41,98.44)	12.86 (10.06,44.26)
OAI-iMorphics-FLASH	Cartilage(Tibia)	104.73 (66.96,174.09)	115.32 (63.55,146.85)	63.73 (47.28,112.59)	97.79 (60.01,102.95)	37.72 (36.13,57.15)	41.05 (35.94,45.42)	70.10 (54.49,84.13)	17.23 (15.22,21.41)	17.23 (15.22,21.41)	19.33 (13.86,31.77)	17.23 (15.65,23.25)
OAI-iMorphics-FLASH	Patella	233.53 (202.77,270.39)	61.31 (10.54,202.30)	8.74 (7.84,10.78)	271.09 (257.78,285.27)	161.86 (110.86,202.35)	36.99 (24.60,37.61)	16.27 (11.62,20.19)	11.20 (9.17,11.41)	23.42 (10.23,140.33)	34.60 (11.38,55.23)	37.26 (10.35,49.76)
OAI-iMorphics-FLASH	Cartilage(Patella)	254.17 (215.31,304.35)	235.96 (180.28,279.54)	175.97 (147.65,207.13)	193.33 (147.65,207.13)	49.29 (34.60,62.93)	54.06 (50.04,60.07)	11.20 (9.44,113.16)	9.42 (6.34,10.59)	10.10 (8.28,15.29)	10.11 (7.25,14.94)	10.02 (8.76,11.96)
OAI-AKOA	Bone(Femur)	37.91 (36.46,39.36)	36.76 (35.41,38.11)	56.79 (47.69,65.90)	71.21 (55.49,86.93)	65.37 (54.15,76.58)	39.55 (35.34,43.76)	42.34 (39.66,45.01)	20.47 (17.20,23.75)	13.57 (10.69,16.46)	17.16 (12.22,10.21)	21.98 (14.63,29.33)
OAI-AKOA	Bone(Tibia)	71.08 (68.35,73.81)	61.68 (61.64,61.71)	63.47 (61.67,65.88)	195.82 (194.86,196.78)	57.96 (54.98,60.94)	35.17 (33.78,36.57)	27.22 (25.84,28.60)	7.60 (7.37,7.83)	36.45 (26.19,46.78)	7.21 (7.14,7.28)	29.70 (21.12,28.99)
OAI-AKOA	Cartilage(Femur)	120.22 (81.76,158.69)	113.07 (97.38,128.76)	65.54 (64.75,66.34)	93.51 (90.81,96.21)	85.67 (84.13,87.20)	83.14 (82.16,84.12)	65.56 (64.46,66.65)	38.82 (25.09,52.55)	18.29 (14.62,21.96)	15.66 (12.33,18.99)	39.88 (26.72,53.04)
OAI-AKOA	Cartilage(Patella)	178.48 (171.34,185.63)	145.51 (138.69,152.33)	162.65 (158.43,166.88)	159.38 (136.39,197.49)	49.17 (44.03,54.31)	59.94 (57.28,62.60)	8.80 (8.43,9.17)	55.08 (53.27,56.89)	35.05 (22.55,47.55)	10.22 (8.64,11.27)	20.55 (15.22,25.97)
OAI-AKOA	Cartilage(Tibia)	94.45 (94.42,94.47)	76.85 (67.24,86.46)	75.74 (66.60,84.87)	37.34 (36.73,37.94)	68.54 (63.76,73.33)	63.61 (49.62,77.60)	32.88 (23.15,42.61)	23.53 (17.60,29.46)	14.73 (13.22,16.24)	16.55 (14.70,18.40)	30.28 (26.49,40.47)
OAI-AKOA	Patella	228.13 (198.56,257.70)	240.20 (237.89,242.51)	58.89 (33.48,84.30)	276.48 (254.50,298.46)	42.77 (36.38,49.15)	26.67 (24.30,29.03)	11.78 (11.39,12.18)	7.61 (7.48,7.74)	6.94 (6.67,7.21)	8.43 (8.40	

Dataset	Target	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2dloss	GVT-v1	GVT-v2	SAMRI-3D
AMOS-MR	Aorta	11.24 (5.36,30.81)	10.16 (4.55,30.04)	18.10 (10.94,28.57)	13.07 (5.55,29.14)	31.73 (16.53,50.80)	26.35 (23.30,47.95)	7.47 (4.69,87.66)	6.85 (3.76,68.66)	10.48 (4.55,35.26)	7.87 (4.49,33.74)	7.69 (3.91,81.87)
AMOS-MR	Duodenum	81.16 (62.91,105.08)	61.19 (49.47,97.50)	76.11 (69.44,112.12)	67.13 (52.10,68.37)	38.68 (35.86,54.71)	27.05 (19.85,29.35)	26.16 (18.19,35.04)	37.46 (33.75,43.66)	52.30 (45.44,65.58)	39.32 (34.54,87.72)	37.30 (33.64,65.58)
AMOS-MR	Esophagus	62.60 (24.31,63.94)	60.25 (33.43,77.04)	59.06 (56.23,80.09)	101.05 (95.54,123.83)	51.98 (38.11,62.82)	16.15 (12.33,19.01)	8.85 (7.41,57.22)	11.02 (6.79,13.53)	21.54 (6.89,63.96)	13.69 (6.13,19.01)	18.30 (8.49,22.99)
AMOS-MR	Gallbladder	47.26 (23.97,146.41)	60.12 (29.16,138.16)	20.41 (17.86,45.01)	75.43 (47.33,102.48)	17.23 (8.46,26.75)	12.39 (11.97,17.79)	7.31 (6.37,8.13)	10.00 (8.34,12.56)	11.09 (9.12,28.34)	7.50 (6.49,9.87)	10.97 (5.54,20.64)
AMOS-MR	Inferior Vena Cava	63.11 (61.35,66.79)	57.47 (50.60,82.63)	75.93 (39.39,93.37)	56.83 (47.28,65.62)	49.55 (43.24,55.27)	20.43 (19.75,20.92)	7.75 (7.38,11.70)	18.27 (8.76,48.36)	16.24 (14.35,20.15)	27.72 (13.48,37.06)	10.25 (9.13,24.66)
AMOS-MR	Left Adrenal Gland	61.43 (35.91,67.94)	81.58 (61.66,105.79)	51.56 (28.80,103.69)	58.53 (49.70,87.85)	22.47 (9.55,37.21)	14.03 (11.72,18.05)	10.44 (8.18,12.66)	5.50 (4.43,41.81)	6.56 (6.10,14.80)	7.53 (4.10,9.22)	5.91 (4.39,13.95)
AMOS-MR	Left Kidney	35.52 (31.16,64.59)	55.80 (39.87,79.07)	146.59 (75.60,172.11)	105.81 (72.95,115.29)	28.13 (18.58,34.83)	5.95 (5.19,16.28)	112.26 (105.22,149.78)	4.99 (3.31,16.90)	11.47 (4.17,22.52)	6.85 (4.48,15.84)	9.70 (6.22,17.83)
AMOS-MR	Liver	34.33 (19.38,66.72)	47.08 (33.22,60.78)	30.92 (26.17,50.13)	69.61 (66.50,72.08)	54.59 (47.69,63.92)	55.03 (51.21,58.32)	73.31 (45.68,90.81)	24.81 (19.61,53.18)	24.81 (19.61,57.83)	55.48 (27.28,61.60)	26.11 (22.86,39.91)
AMOS-MR	Pancreas	64.28 (51.99,79.21)	56.08 (35.10,87.29)	51.69 (46.19,80.20)	78.40 (38.96,121.46)	43.50 (37.39,48.93)	32.50 (22.93,56.03)	40.32 (26.44,78.92)	18.86 (15.87,27.14)	29.85 (14.70,39.93)	20.54 (17.93,32.55)	17.48 (13.50,84.43)
AMOS-MR	Right Adrenal Gland	78.42 (46.93,95.41)	74.03 (32.00,132.87)	87.91 (50.92,95.26)	93.39 (70.07,105.58)	10.00 (7.94,28.58)	19.99 (12.77,25.89)	9.28 (6.32,11.17)	17.08 (12.63,24.96)	15.96 (11.26,29.62)	7.63 (5.61,10.32)	8.70 (6.61,9.95)
AMOS-MR	Right Kidney	46.13 (29.05,68.63)	48.57 (29.51,68.08)	42.21 (19.93,81.16)	33.98 (16.47,87.39)	36.25 (21.91,59.87)	8.13 (6.44,25.96)	12.33 (8.93,77.58)	5.74 (4.52,9.08)	14.38 (7.04,21.84)	4.47 (3.90,6.01)	5.74 (3.95,10.96)
AMOS-MR	Spleen	60.74 (34.98,80.05)	117.99 (68.86,146.50)	63.55 (55.71,76.04)	99.86 (73.49,123.64)	39.40 (31.45,48.94)	18.95 (14.31,25.29)	60.99 (21.62,101.31)	7.50 (3.87,19.07)	6.95 (3.58,24.66)	29.69 (12.89,45.25)	23.33 (12.94,31.63)
AMOS-MR	Stomach	38.73 (29.93,68.54)	51.95 (38.40,61.48)	73.40 (60.44,106.00)	62.78 (45.90,101.82)	51.77 (36.93,91.19)	39.18 (33.19,53.39)	19.19 (13.97,70.59)	25.66 (18.37,35.69)	38.95 (34.35,40.82)	35.43 (19.13,58.88)	48.07 (39.78,60.74)
MSK-Shoulder	Humerus cartilage	60.44 (31.14,120.02)	129.25 (127.84,159.67)	95.36 (28.30,108.30)	159.54 (135.78,203.26)	48.18 (43.84,54.05)	57.36 (52.50,63.03)	17.58 (15.00,28.53)	25.04 (20.10,46.70)	53.24 (24.76,106.12)	13.60 (10.82,14.76)	14.76 (14.56,15.17)
MSK-Shoulder	Scapula cartilage	110.15 (105.55,130.50)	127.74 (106.81,220.44)	163.80 (126.94,290.24)	156.06 (102.84,157.11)	80.02 (61.23,81.00)	94.47 (88.82,96.34)	82.02 (70.22,88.28)	119.12 (63.21,123.50)	60.96 (37.75,75.17)	107.05 (76.82,142.56)	111.95 (55.57,134.56)
CervicalCancer	Cervix Cancer	17.61 (17.61,17.61)	15.81 (15.81,15.81)	12.53 (12.53,12.53)	16.25 (16.25,16.25)	17.12 (17.12,17.12)	13.04 (13.04,13.04)	14.87 (14.87,14.87)	16.19 (16.19,16.19)	15.84 (15.84,15.84)	15.84 (15.84,15.84)	16.79 (16.79,16.79)
Mean $\pm$ SD		67.73 $\pm$ 56.64	66.33 $\pm$ 48.99	58.82 $\pm$ 45.55	87.03 $\pm$ 60.69	45.84 $\pm$ 30.03	35.38 $\pm$ 20.58	32.94 $\pm$ 28.99	20.93 $\pm$ 22.76	22.96 $\pm$ 22.03	21.41 $\pm$ 22.70	23.53 $\pm$ 25.51

Table S11. MSD by target — median (Q1, Q3) across volumes; final row is mean $\pm$ SD.

Dataset	Target	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2dloss	GVT-v1	GVT-v2	SAMRI-3D
MSD-BrainTumour	Edema	9.46 (7.59,12.05)	8.77 (6.69,11.26)	8.76 (6.51,11.72)	9.74 (6.94,12.21)	9.68 (7.55,12.42)	3.77 (2.82,5.24)	3.16 (2.56,3.63)	6.68 (3.95,9.03)	4.83 (3.10,7.19)	6.20 (3.62,9.10)	5.71 (3.65,7.94)
MSD-BrainTumour	Enhancing tumour	13.54 (10.43,16.43)	11.76 (7.16,15.23)	10.36 (6.63,15.21)	12.61 (8.50,17.21)	7.44 (5.17,12.06)	4.44 (3.09,6.46)	2.22 (1.87,2.68)	4.94 (3.30,7.76)	3.85 (2.60,5.51)	5.46 (3.30,7.68)	4.79 (3.07,7.03)
MSD-BrainTumour	Non-enhancing tumour	14.58 (10.72,20.55)	11.91 (7.75,16.86)	10.51 (6.92,15.60)	13.91 (9.38,18.56)	8.87 (5.92,13.71)	4.44 (3.20,6.04)	2.27 (1.97,2.84)	7.17 (4.78,10.30)	6.46 (4.64,9.58)	7.17 (4.99,10.75)	6.88 (4.73,10.53)
BraTS-FLAIR	Brain Tumour	3.32 (2.54,4.55)	3.30 (2.58,4.97)	3.18 (2.23,4.36)	3.26 (2.20,4.94)	4.73 (3.54,6.01)	2.36 (2.05,2.79)	2.42 (2.08,2.91)	2.24 (1.74,2.94)	2.24 (1.73,3.02)	1.93 (1.60,2.64)	2.07 (1.65,3.07)
BraTS-T1	Brain Tumour	4.45 (3.52,5.82)	4.64 (3.73,6.01)	4.83 (3.60,6.04)	5.23 (3.81,7.34)	5.50 (4.44,7.11)	4.13 (3.57,5.23)	3.02 (2.60,3.45)	3.29 (2.60,4.29)	3.46 (2.58,4.32)	2.88 (2.46,3.75)	3.31 (2.49,4.02)
BraTS-T1CE	Brain Tumour	2.32 (1.68,3.18)	2.62 (1.68,3.27)	2.31 (1.69,3.31)	2.37 (1.78,3.29)	2.48 (1.94,3.13)	3.14 (2.70,3.73)	1.90 (1.71,2.22)	1.90 (1.71,2.22)	1.55 (1.26,2.15)	1.47 (1.25,2.40)	1.67 (1.32,2.74)
ISLES2022-DWI	Ischemic Stroke	3.91 (1.63,6.38)	3.56 (2.20,6.30)	3.22 (1.43,5.91)	4.95 (1.86,6.41)	1.53 (1.15,3.21)	6.25 (4.78,7.63)	1.17 (1.10,1.30)	1.29 (1.05,1.83)	1.32 (1.07,1.82)	1.18 (1.03,1.67)	1.29 (1.08,1.64)
ISLES2022-ADC	Ischemic Stroke	3.25 (1.49,5.49)	4.26 (1.99,6.80)	2.92 (1.57,4.48)	3.85 (2.05,6.62)	2.82 (1.80,3.64)	5.28 (4.13,6.66)	1.22 (1.16,1.42)	1.49 (1.17,2.33)	1.60 (1.18,2.22)	1.50 (1.19,2.52)	1.42 (1.17,2.24)
WMH-T1	White Matter Hyperintensity	9.01 (5.90,12.37)	8.32 (5.70,12.45)	6.99 (3.70,12.11)	9.91 (6.22,15.05)	2.50 (1.35,5.60)	6.93 (6.01,8.11)	1.68 (1.38,3.04)	1.85 (1.31,2.68)	2.06 (1.34,3.30)	1.51 (1.21,2.15)	1.54 (1.22,2.43)
WMH-FLAIR	White Matter Hyperintensity	10.19 (5.87,20.96)	9.00 (5.39,13.84)	7.38 (3.63,15.91)	10.20 (6.08,16.12)	2.36 (1.33,4.58)	6.87 (5.71,8.03)	1.76 (1.36,2.81)	1.83 (1.27,3.89)	3.07 (1.48,6.21)	1.71 (1.18,3.65)	1.61 (1.21,3.52)
PI-CAI	Prostate	3.42 (2.36,4.74)	3.79 (3.09,5.32)	4.08 (3.09,6.00)	4.84 (3.75,7.74)	4.61 (3.28,5.70)	9.24 (7.78,11.02)	2.04 (1.74,2.19)	1.49 (1.29,1.64)	1.68 (1.40,2.01)	1.58 (1.46,1.89)	1.58 (1.41,1.79)
Prostate-T2	Prostate	2.09 (1.75,2.50)	2.43 (1.93,2.94)	2.39 (1.91,2.75)	2.37 (2.12,2.88)	3.35 (2.65,5.06)	9.19 (7.16,11.87)	1.47 (1.38,1.73)	1.48 (1.33,1.86)	1.48 (1.38,1.75)	1.50 (1.38,1.62)	1.50 (1.38,1.62)
Prostate-ADC	Prostate	1.66 (1.39,1.95)	1.56 (1.32,1.83)	1.58 (1.42,1.75)	1.50 (1.32,1.72)	1.94 (1.56,2.52)	2.21 (1.96,2.42)	1.29 (1.19,1.55)	1.12 (1.08,1.28)	1.14 (1.07,1.25)	1.15 (1.08,1.32)	1.14 (1.10,1.25)
QIN-Prostate	Prostate	3.42 (3.01,3.66)	2.51 (2.04,3.57)	3.86 (2.85,5.11)	3.02 (2.53,3.66)	3.26 (2.57,4.21)	13.08 (5.46,17.60)	2.37 (2.29,3.86)	1.58 (1.53,2.50)	1.63 (1.57,2.75)	1.91 (1.59,2.80)	1.64 (1.51,2.43)
PROMISE12	Prostate	4.97 (3.46,7.85)	5.82 (4.96,8.29)	8.83 (6.66,12.46)	8.22 (7.07,10.64)	4.77 (3.29,6.78)	6.55 (5.49,10.95)	1.72 (1.58,1.88)	1.45 (1.34,1.52)	1.80 (1.58,2.19)	1.45 (1.33,1.66)	1.47 (1.40,1.78)
MSD-Prostate	Prostate Central Zone	5.18 (4.62,5.82)	5.04 (4.02,6.19)	5.62 (3.71,6.45)	4.53 (3.87,5.09)	4.81 (3.59,6.75)	3.54 (3.45,3.85)	1.54 (1.45,1.67)	1.95 (1.85,2.06)	2.17 (2.02,2.31)	2.03 (1.93,2.23)	2.43 (2.03,2.75)
MSD-Prostate	Prostate Peripheral Zone	5.64 (3.58,6.56)	5.12 (3.65,5.57)	5.68 (3.80,8.28)	5.41 (3.90,6.11)	5.28 (4.31,7.59)	5.22 (4.64,5.98)	2.66 (2.47,2.75)	2.77 (2.03,3.64)	2.90 (2.37,3.18)	2.63 (2.15,3.90)	2.43 (2.03,2.75)
QIN-Lesion	Prostate Peripheral Zone	2.13 (1.43,3.19)	2.17 (1.41,2.44)	2.63 (2.03,2.99)	2.11 (1.62,2.44)	2.24 (1.56,2.47)	6.14 (4.90,6.60)	2.77 (1.61,3.74)	1.53 (1.41,1.84)	1.53 (1.37,1.84)	1.74 (1.34,2.19)	1.74 (1.34,2.19)
ZIB-OAI	Bone(Femur)	4.03 (2.30,10.10)	3.44 (2.02,4.67)	2.38 (1.94,8.59)	4.97 (2.40,8.60)	15.41 (9.03,24.67)	7.72 (6.46,9.36)	7.65 (7.19,8.00)	1.57 (1.42,1.87)	1.56 (1.42,1.77)	1.55 (1.41,1.75)	1.80 (1.53,2.79)
ZIB-OAI	Bone(Tibia)	2.50 (1.98,8.48)	6.59 (5.80,8.13)	6.47 (4.97,7.83)	7.41 (4.79,12.11)	17.36 (15.04,24.77)	7.88 (6.13,8.65)	7.37 (1.26,5.12)	4.59 (1.48,5.43)	1.42 (1.33,1.82)	1.64 (1.37,2.58)	1.64 (1.37,2.58)
ZIB-OAI	Cartilage(Femur)	11.43 (8.75,14.71)	10.23 (8.38,13.13)	12.03 (9.79,14.99)	12.11 (9.85,15.60)	22.23 (20.23,27.23)	18.78 (17.88,20.30)	12.47 (11.28,14.26)	1.73 (1.49,2.38)	1.65 (1.46,1.99)	1.57 (1.40,1.93)	1.76 (1.54,2.21)
ZIB-OAI	Cartilage(Tibia)	9.16 (5.88,13.07)	11.96 (9.15,15.94)	12.08 (9.92,17.70)	13.55 (9.93,17.73)	18.11 (8.37,20.52)	18.11 (8.37,20.52)	2.90 (2.70,3.06)	2.73 (2.71,4.31)	2.81 (2.73,4.31)	2.52 (2.24,2.97)	2.91 (2.59,3.49)
OAI-iMorphics-DESS	Bone(Femur)	4.09 (2.91,12.27)	3.18 (2.12,4.96)	3.73 (2.13,10.24)	4.48 (2.49,7.58)	11.98 (7.15,25.43)	12.07 (8.32,14.54)	7.33 (6.94,7.66)	1.60 (1.50,1.87)	1.55 (1.42,1.85)	1.56 (1.49,1.69)	1.81 (1.64,2.43)
OAI-iMorphics-DESS	Cartilage(Femur)	9.65 (7.08,13.97)	9.63 (8.64,11.50)	11.42 (9.52,14.72)	10.91 (8.87,14.78)	24.99 (19.85,28.16)	19.63 (17.72,22.25)	13.83 (13.35,14.54)	1.96 (1.47,3.07)	1.60 (1.42,2.07)	1.62 (1.39,2.11)	1.77 (1.59,2.35)
OAI-iMorphics-DESS	Bone(Tibia)	2.36 (2.03,6.72)	5.73 (2.28,7.72)	6.59 (3.59,7.75)	7.07 (4.61,8.68)	17.75 (12.83,25.76)	15.41 (12.64,18.29)	6.75 (6.34,7.62)	1.29 (1.24,1.38)	1.76 (1.27,5.24)	1.30 (1.25,1.44)	1.39 (1.26,1.85)
OAI-iMorphics-DESS	Cartilage(Tibia)	10.69 (7.01,19.09)	15.20 (10.92,20.48)	16.66 (10.90,23.89)	12.01 (7.69,15.30)	18.00 (13.93,20.48)	8.20 (5.93,10.15)	2.75 (2.50,2.86)	2.09 (2.03,2.30)	2.20 (2.38,2.95)	2.49 (2.17,2.98)	2.97 (2.40,3.74)
OAI-iMorphics-DESS	Patella	50.95 (38.69,59.18)	23.40 (1.66,39.39)	1.41 (1.31,1.66)	41.40 (5.11,49.21)	14.97 (12.82,18.59)	5.97 (4.99,7.50)	3.41 (3.22,3.87)	1.36 (1.28,1.64)	1.49 (1.28,9.26)	1.46 (1.28,2.19)	3.52 (1.35,9.05)
OAI-iMorphics-DESS	Cartilage(Patella)	38.43 (34.51,47.99)	28.47 (19.03,37.92)	25.61 (6.69,33.73)	32.30 (20.07,43.18)	11.23 (9.97,13.29)	10.88 (8.99,11.99)	2.47 (2.23,3.83)	1.48 (1.33,1.83)	1.48 (1.30,1.74)	1.32 (1.23,1.63)	1.34 (1.21,1.63)
OAI-iMorphics-TSE	Bone(Femur)	9.60 (6.81,10.06)	5.15 (3.72,7.80)	8.68 (7.12,9.86)	8.73 (7.46,10.07)	5.39 (4.25,8.01)	8.53 (7.11,10.66)	6.51 (6.05,7.09)	2.97 (2.55,5.79)	3.00 (2.54,5.13)	3.06 (2.34,5.13)	5.17 (2.63,6.28)
OAI-iMorphics-TSE	Cartilage(Femur)	15.10 (13.10,16.20)	10.01 (5.12,12.33)	16.93 (14.43,20.30)	14.96 (12.65,16.03)	18.43 (16.66,23.44)	14.48 (12.29,16.12)	12.35 (11.77,12.83)	8.41 (3.19,14.35)	4.07 (2.26,11.96)	5.65 (2.39,11.76)	9.32 (4.40,13.63)
OAI-iMorphics-TSE	Bone(Tibia)	3.31 (2.95,4.34)	5.13 (4.05,7.44)	4.87 (3.38,8.66)	17.67 (11.88,25.78)	9.41 (6.68,11.55)	6.73 (5.82,7.30)	5.84 (4.79,6.23)	2.85 (2.43,3.16)	3.30 (2.99,4.09)	2.79 (2.60,3.19)	3.04 (2.64,3.60)
OAI-iMorphics-TSE	Cartilage(Tibia)	5.06 (4.04,8.42)	5.61 (4.47,9.67)	7.75 (6.31,10.08)	7.68 (6.29,12.78)	8.66 (7.82,10.15)	8.66 (7.82,10.15)	2.80 (2.61,2.94)	3.20 (2.63,8.02)	2.76 (2.40,3.54)	3.01 (2.27,6.02)	3.30 (2.42,4.14)
OAI-iMorphics-TSE	Patella	22.20 (7.96,32.60)	4.90 (2.85,13.68)	3.48 (2.40,4.59)	35.30 (6.11,41.99)	6.41 (5.09,9.09)	6.97 (4.90,8.54)	3.27 (2.96,3.61)	2.40 (2.14,3.03)	2.59 (2.24,3.17)	2.27 (1.85,3.14)	2.27 (1.86,2.85)
OAI-iMorphics-TSE	Cartilage(Patella)	13.15 (8.95,17.68)	6.04 (3.95,17.68)	14.51 (10.59,18.90)	18.95 (13.95,27.31)	9.46 (6.45,11.46)	10.23 (8.79,11.20)	2.73 (2.45,3.20)	2.02 (1.82,3.20)	2.31 (1.82,5.25)	2.17 (1.86,3.60)	2.24 (1.90,2.76)
OAI-iMorphics-FLASH	Bone(Femur)	13.08 (5.21,15.12)	6.04 (3.71,8.98)	2.03 (1.65,9.45)	8.11 (7.37,10.93)	3.25 (2.52,4.47)	14.42 (13.03,20.15)	6.85 (6.70,7.22)	1.40 (1.32,1.61)	1.65 (1.44,4.42)	1.73 (1.40,5.43)	2.09 (1.46,5.05)
OAI-iMorphics-FLASH	Cartilage(Femur)	12.95 (11.61,18.23)	12.31 (8.45,15.22)	10.51 (8.58,14.82)	9.65 (7.48,11.97)	17.04 (16.65,19.57)	19.78 (16.91,22.74)	12.51 (11.83,12.89)	2.59 (1.99,6.61)	1.66 (1.59,1.81)	1.80 (1.73,1.95)	2.00 (1.76,5.57)
OAI-iMorphics-FLASH	Bone(Tibia)	31.30 (15.57,45.67)	6.07 (5.61,7.33)	8.12 (3.29,32.94)	13.33 (11.76,32.14)	5.33 (4.44,8.87)	9.34 (7.74,9.93)	5.45 (5.11,6.10)	1.57 (1.32,1.75)	5.82 (4.94,34.98)	2.70 (1.43,10.71)	1.65 (1.29,2.19)
OAI-iMorphics-FLASH	Cartilage(Tibia)	9.47 (3.87,15.90)	6.75 (4.80,10.31)	7.34 (4.68,11.88)	7.09 (6.12,12.41)	3.99 (3.38,10.67)	7.83 (5.75,8.55)	3.43 (3.20,3.66)	2.22 (2.04,2.55)	2.46 (2.24,3.54)	2.73 (2.28,3.74)	2.81 (2.50,3.81)
OAI-iMorphics-FLASH	Patella	39.80 (34.03,53.18)	13.03 (1.57,34.00)	1.48 (1.40,1.54)	53.03 (41.83,59.51)	20.77 (17.50,35.75)	4.45 (3.57,5.59)	3.24 (2.91,3.59)	1.50 (1.42,1.72)	2.67 (1.46,16.24)	3.81 (1.66,6.11)	4.48 (1.45,6.18)
OAI-iMorphics-FLASH	Cartilage(Patella)	46.55 (41.10,56.63)	38.51 (28.16,47.67)	26.08 (12.26,41.83)	31.54 (10.74,42.76)	7.46 (3.36,8.25)	6.53 (5.66,9.34)	1.99 (1.59,2.34)	1.86 (1.53,2.41)	1.46 (1.44,1.86)	1.60 (1.44,2.12)	1.60 (1.44,2.12)
OAI-AKOA	Bone(Femur)	5.08 (4.63,5.54)	4.34 (4.20,4.48)	10.50 (8.45,12.56)	7.78 (7.27,8.28)	7.47 (6.47,8.47)	7.60 (7.21,7.98)	6.10 (6.08,6.11)	2.32 (1.83,2.80)	2.05 (1.74,2.35)	2.27 (1.80,2.73)	2.79 (2.02,3.57)
OAI-AKOA	Bone(Tibia)	2.01 (1.97,2.04)	8.01 (7.88,8.14)	6.19 (5.71,6.68)	12.43 (11.05,13.82)	12.09 (11.77,12.41)	5.49 (5.37,5.60)	5.75 (5.58,5.91)	1.32 (1.32,1.32)	3.48 (2.42,5.44)	1.31 (1.31,1.32)	1.57 (1.54,1.60)
OAI-AKOA	Cartilage(Femur)	10.08 (6.64,13.52)	10.48 (10.27,10.69)	9.15 (6.95,11.35)	13.72 (12.58,14.86)	15.99 (15.46,16.52)	17.18 (16.53,17.84)	12.15 (12.11,12.19)	5.00 (3.14,6.85)	2.21 (1.75,6.27)	1.66 (1.46,1.86)	5.25 (3.37,7.13)
OAI-AKOA	Cartilage(Patella)	34.14 (34.13,34.15)	27.86 (25.50,30.23)	30.02 (27.47,32.58)	36.19 (32.17,19.12)	12.59 (10.51,14.68)	12.25 (11.98,12.52)	1.90 (1.81,1.99)	4.40 (3.49,4.41)	4.16 (2.78,5.55)	1.46 (1.45,1.47)	2.81 (2.13,5.50)
OAI-AKOA	Cartilage(Tibia)	9.49 (9.42,9.57)	10.53 (8.90,12.16)	11.37 (10.14,12.61)	3.82 (3.82,3.83)	16.03 (14.38,17.67)	11.33 (9.76,12.91)	2.63 (2.62,2.64)	3.04 (2.43,3.64)	2.32 (2.03,6.21)	2.40 (2.13,2.66)	3.33 (2.88,3.79)
OAI-AKOA	Patella	50.73 (46.55,54.90)	55.75 (50.87,60.63)	7.01 (7.66,20.37)	63.31 (61.44,65.18)	8.78 (7.91,9.64)	2.78 (2.62,2.94)	7.91 (7.25,8.44)	1.28 (1.26,1.30)	1.26 (1.25,1.27)	1.38 (1.36,1.39)	1.37 (1.35,1.39)
MSK-T2	Bone(Femur)	2.35 (2.27,2.75)	2.41 (2.24,3.07)	2.36 (2.31,2.51)	4.80 (2.30,7.78)	9.29 (8.57,10.77)	8.77 (7.99,10.48)	7.25 (7.15,7.64)	1.86 (1.85,2.01)	1.85 (1.80,2.06)	1.86 (1.82,2.76)	2.26 (1.90,3.00)
MSK-T2	Bone(Tibia)	6.97 (1.82,13.26)	10.49 (3.53,17.16)	1.83 (1.77,1.86)	9.78 (1.82,19.88)	22.96 (15.81,28.91)	4.25 (3.47,5.07)	6.24 (5.71,6.44)	3.38 (3.68,6.37)	8.70 (2.75,14.51)	13.50 (10.74,14.49)	10.10 (5.93,13.90)
MSK-T2	Cartilage(Femur)	6.55 (6.06,6.86)	6.76 (5.81,7.77)	4.01 (3.67,6.43)	5.11 (4.11,5.97)	15.76 (12.39,24.84)	15.67 (15.09,18.39)	9.42 (9.30,9.59)	1.24 (1.16,1.39)	1.34 (1.17,1.53)	1.19 (1.17,1.21)	1.56 (1.21,1.93)
MSK-T2	Cartilage(Patella)	1.34 (1.25,1.43)	1.80 (1.37,2.69)	1.85 (1.44,2.24)	15.32 (9.10,34.06)	8.62 (6.32,11.34)	7.45 (6.80,8.14)	1.95 (1.82,2.15)	1.31 (1.28,1.33)	1.41 (1.33,1.44)	1.28 (1.20,1.37)	1.28 (1.19,1.37)
MSK-T2	Cartilage(Tibia)	17.33 (15.47,18.76)	16.68 (14.27,19.53)	13.88 (10.84,16.21)	17.55 (16.58,19.65)	5.29 (4.97,6.47)	5.71 (4.92,7.01)	2.92 (2				

Dataset	Target	SAM2-Tiny	SAM2-Small	SAM2-Base+	SAM2-Large	Medical-SAM2	SAM-Med3D	SAMed-2	SAMRI-3D-2dloss	GVT-v1	GVT-v2	SAMRI-3D
AMOS-MR	Aorta	1.41 (1.19,4.72)	1.43 (1.23,4.76)	1.52 (1.23,3.75)	1.75 (1.26,5.37)	5.63 (2.74,15.00)	6.71 (5.45,7.62)	1.39 (1.34,1.66)	1.37 (1.19,12.67)	1.34 (1.19,3.99)	1.33 (1.16,4.96)	1.34 (1.16,15.17)
AMOS-MR	Duodenum	15.89 (11.78,21.55)	11.42 (10.26,22.35)	13.06 (10.84,17.79)	13.93 (11.29,16.85)	8.99 (7.12,11.15)	7.05 (6.67,7.76)	3.72 (2.73,5.03)	7.27 (5.34,8.60)	9.06 (8.49,9.76)	9.74 (6.08,13.83)	7.11 (4.59,9.42)
AMOS-MR	Esophagus	10.92 (4.18,15.26)	10.30 (6.99,16.03)	13.27 (11.52,26.43)	18.07 (14.62,28.81)	10.64 (9.55,10.85)	3.96 (3.70,4.36)	1.58 (1.45,2.38)	2.59 (1.72,2.96)	4.95 (1.71,8.23)	1.91 (1.54,3.28)	2.86 (1.72,3.58)
AMOS-MR	Gallbladder	6.93 (3.09,10.51)	13.46 (6.93,16.73)	5.90 (2.97,15.96)	15.36 (6.45,17.23)	3.00 (1.60,4.25)	3.86 (3.22,5.88)	1.64 (1.51,1.88)	3.20 (2.31,3.51)	3.35 (2.66,9.01)	2.62 (2.16,3.10)	3.16 (1.77,6.04)
AMOS-MR	Inferior Vena Cava	13.37 (12.74,19.17)	8.55 (5.35,16.46)	16.28 (9.86,20.32)	14.03 (5.93,21.68)	11.31 (7.86,12.67)	5.61 (5.14,6.20)	1.67 (1.49,1.89)	3.15 (1.77,10.32)	2.62 (2.15,2.69)	4.14 (1.83,5.90)	1.65 (1.58,3.99)
AMOS-MR	Left Adrenal Gland	9.09 (5.72,16.77)	11.31 (10.04,14.15)	11.29 (7.16,19.81)	13.02 (10.83,17.16)	4.59 (2.64,9.67)	3.81 (3.52,4.45)	1.86 (1.72,2.04)	1.63 (1.34,9.18)	1.63 (1.44,2.39)	1.85 (1.43,2.20)	1.60 (1.29,2.93)
AMOS-MR	Left Kidney	6.23 (4.37,11.22)	6.45 (5.99,8.32)	9.28 (5.00,11.55)	6.47 (4.33,8.85)	5.15 (3.67,7.01)	2.54 (1.82,3.59)	2.41 (2.30,2.63)	1.45 (1.20,2.40)	2.30 (1.43,2.61)	1.78 (1.47,3.85)	2.28 (1.46,3.14)
AMOS-MR	Liver	3.07 (1.84,5.37)	4.87 (4.23,6.30)	3.60 (3.08,4.12)	4.87 (3.86,6.63)	7.03 (4.83,9.97)	13.99 (12.47,17.83)	5.41 (4.57,7.87)	2.53 (1.60,3.98)	2.65 (1.63,4.65)	4.12 (2.20,4.98)	2.59 (1.82,3.20)
AMOS-MR	Pancreas	12.20 (9.72,15.91)	12.56 (8.08,15.91)	12.48 (8.53,17.57)	12.20 (8.17,16.27)	10.86 (8.92,13.35)	10.04 (5.97,13.26)	4.07 (3.53,5.59)	2.98 (2.27,5.34)	4.46 (2.01,7.29)	2.96 (2.36,3.87)	2.79 (2.16,13.12)
AMOS-MR	Right Adrenal Gland	17.79 (11.88,22.59)	15.05 (8.22,26.41)	18.71 (11.82,24.17)	17.86 (16.72,21.52)	2.12 (1.77,3.10)	5.03 (4.29,6.16)	1.96 (1.70,2.13)	3.16 (2.27,3.79)	2.79 (2.39,4.56)	1.77 (1.53,2.39)	1.86 (1.62,2.51)
AMOS-MR	Right Kidney	5.70 (4.16,7.39)	6.65 (4.33,9.05)	7.52 (4.54,16.65)	6.14 (3.34,7.81)	6.47 (4.48,12.89)	2.18 (2.00,4.57)	2.26 (2.08,2.44)	2.00 (1.78,2.22)	4.00 (2.57,5.03)	1.91 (1.68,2.18)	2.19 (1.70,2.21)
AMOS-MR	Spleen	4.54 (4.34,9.71)	6.08 (4.04,6.87)	8.24 (3.54,11.51)	9.71 (6.25,12.45)	5.53 (5.08,8.37)	4.23 (2.40,5.57)	3.02 (2.76,3.46)	1.22 (1.08,1.35)	1.41 (1.10,1.76)	2.29 (1.28,3.25)	1.49 (1.27,2.21)
AMOS-MR	Stomach	8.62 (3.46,10.47)	7.61 (5.97,7.92)	8.67 (6.04,14.72)	8.88 (6.56,9.96)	13.73 (9.56,14.47)	11.80 (11.65,15.52)	2.47 (2.20,2.63)	3.44 (2.49,4.19)	3.38 (3.23,3.81)	3.87 (2.45,4.92)	5.83 (5.49,8.93)
MSK-Shoulder	Humerus cartilage	10.61 (4.48,16.78)	12.12 (11.49,25.15)	14.46 (4.41,17.30)	20.53 (14.53,32.91)	5.57 (5.01,6.24)	5.43 (4.82,9.67)	2.80 (2.75,3.19)	2.74 (2.19,4.90)	4.97 (3.05,8.26)	2.48 (1.89,2.50)	2.41 (1.96,2.47)
MSK-Shoulder	Scapula cartilage	19.42 (18.88,22.09)	16.05 (11.99,19.99)	33.41 (28.90,48.05)	18.91 (6.24,28.17)	16.03 (12.99,24.61)	16.44 (12.81,16.88)	11.49 (10.57,12.91)	14.85 (11.69,21.02)	5.01 (3.91,11.01)	18.44 (14.79,27.21)	20.48 (10.99,22.15)
CervicalCancer	Cervix Cancer	2.29 (2.29,2.29)	2.39 (2.39,2.39)	2.57 (2.57,2.57)	2.40 (2.40,2.40)	1.99 (1.99,1.99)	3.28 (3.28,3.28)	2.47 (2.47,2.47)	2.12 (2.12,2.12)	2.08 (2.08,2.08)	1.98 (1.98,1.98)	2.06 (2.06,2.06)
Mean $\pm$ SD		10.31 $\pm$ 10.84	9.23 $\pm$ 8.79	8.88 $\pm$ 8.17	12.38 $\pm$ 11.11	8.70 $\pm$ 6.28	7.47 $\pm$ 4.25	4.05 $\pm$ 3.08	3.13 $\pm$ 3.34	3.41 $\pm$ 4.29	3.12 $\pm$ 3.78	3.15 $\pm$ 3.04

H. Statistical Significance

Table S12. Two-sided paired Wilcoxon signed-rank tests with SAMRI-3D as reference; rows grouped by metric and aggregation level.

Comparison (vs SAMRI-3D)	n	median (ours)	median (other)	median Δ	win-rate (ours)	W	p-value
DSC — Sample-level (n=3376)							
SAMRI-3D-2dloss	3376	0.746	0.734	+0.0005	51.7%	2391185	9.75e-14
GVT-v1	3376	0.746	0.740	+0.0004	50.8%	2510184	3.42e-07
GVT-v2	3376	0.746	0.749	+0.0000	49.9%	2747288	2.17e-01
SAMed-2	3376	0.746	0.642	+0.0699	67.7%	1841746	6.47e-71
SAM2-Small	3376	0.746	0.523	+0.1024	89.0%	306357	0.00e+00
HD — Sample-level (n=3376)							
SAMRI-3D-2dloss	3358	16.093	15.000	+0.0000	43.8%	2210750	1.61e-02
GVT-v1	3344	16.078	16.340	+0.0000	49.1%	2247414	6.12e-03
GVT-v2	3362	16.093	15.000	+0.0000	44.2%	2351564	5.85e-01
SAMed-2	3344	16.125	17.916	-0.9698	54.8%	2323244	4.41e-15
SAM2-Small	3363	16.093	49.020	-18.7418	85.7%	404732	0.00e+00
MSD — Sample-level (n=3376)							
SAMRI-3D-2dloss	3333	2.382	2.363	-0.0055	51.4%	2603225	3.53e-03
GVT-v1	3319	2.371	2.497	-0.0071	51.4%	2564690	1.07e-03
GVT-v2	3334	2.383	2.236	+0.0208	44.9%	2505653	2.18e-06
SAMed-2	3300	2.407	2.674	-0.1205	54.5%	2678516	4.13e-01
SAM2-Small	3340	2.382	6.402	-2.4659	87.6%	425281	0.00e+00
DSC — Dataset-level (n=34)							
SAMRI-3D-2dloss	34	0.822	0.806	+0.0005	50.0%	224	2.15e-01
GVT-v1	34	0.822	0.779	+0.0087	67.6%	151	1.12e-02
GVT-v2	34	0.822	0.817	+0.0005	52.9%	266	6.01e-01
SAMed-2	34	0.822	0.632	+0.0907	85.3%	69	2.87e-05
SAM2-Small	34	0.822	0.641	+0.1353	100.0%	0	1.16e-10
HD — Dataset-level (n=34)							
SAMRI-3D-2dloss	31	13.332	16.186	+0.1842	48.4%	225	6.64e-01
GVT-v1	31	13.332	18.746	-0.4063	64.5%	154	6.64e-02
GVT-v2	31	13.332	15.179	-0.3970	58.1%	238	8.54e-01
SAMed-2	29	16.793	21.592	-4.9570	75.9%	70	8.88e-04
SAM2-Small	31	13.332	52.230	-29.2528	96.8%	1	1.86e-09
MSD — Dataset-level (n=34)							
SAMRI-3D-2dloss	27	2.587	2.201	-0.0024	55.6%	176	7.68e-01
GVT-v1	27	2.587	2.688	-0.1288	74.1%	94	2.14e-02
GVT-v2	27	2.587	2.166	+0.0756	40.7%	154	4.13e-01
SAMed-2	26	2.591	2.693	-0.2533	65.4%	110	9.93e-02
SAM2-Small	27	2.587	5.853	-2.7308	100.0%	0	1.49e-08
DSC — Label-level (n=49)							
SAMRI-3D-2dloss	49	0.730	0.739	+0.0090	71.4%	328	4.08e-03
GVT-v1	49	0.730	0.714	+0.0078	69.4%	326	3.81e-03
GVT-v2	49	0.730	0.735	+0.0017	55.1%	523	3.79e-01
SAMed-2	49	0.730	0.680	+0.0695	69.4%	304	1.75e-03
SAM2-Small	49	0.730	0.513	+0.1907	95.9%	9	1.17e-13
HD — Label-level (n=49)							
SAMRI-3D-2dloss	46	17.013	17.051	-0.1571	52.2%	480	5.16e-01
GVT-v1	46	17.013	22.246	-0.4279	58.7%	471	4.55e-01
GVT-v2	46	17.013	19.117	-0.3249	56.5%	504	6.97e-01
SAMed-2	45	17.232	30.377	-5.9215	64.4%	273	5.11e-03
SAM2-Small	46	17.013	56.254	-34.6531	95.7%	21	1.27e-11
MSD — Label-level (n=49)							
SAMRI-3D-2dloss	44	2.732	2.849	-0.0323	54.5%	440	5.28e-01
GVT-v1	44	2.732	3.382	-0.1768	68.2%	320	4.09e-02
GVT-v2	44	2.732	2.844	-0.0352	52.3%	477	8.40e-01
SAMed-2	43	2.787	2.850	-0.0389	51.2%	465	9.29e-01
SAM2-Small	44	2.732	8.024	-4.7615	95.5%	33	3.65e-10

I. Additional Qualitative Examples

Representative successful (Fig. S2) and failure (Fig. S3) cases across anatomical domains. Ground truth is shown in yellow and the SAMRI-3D prediction in pink.

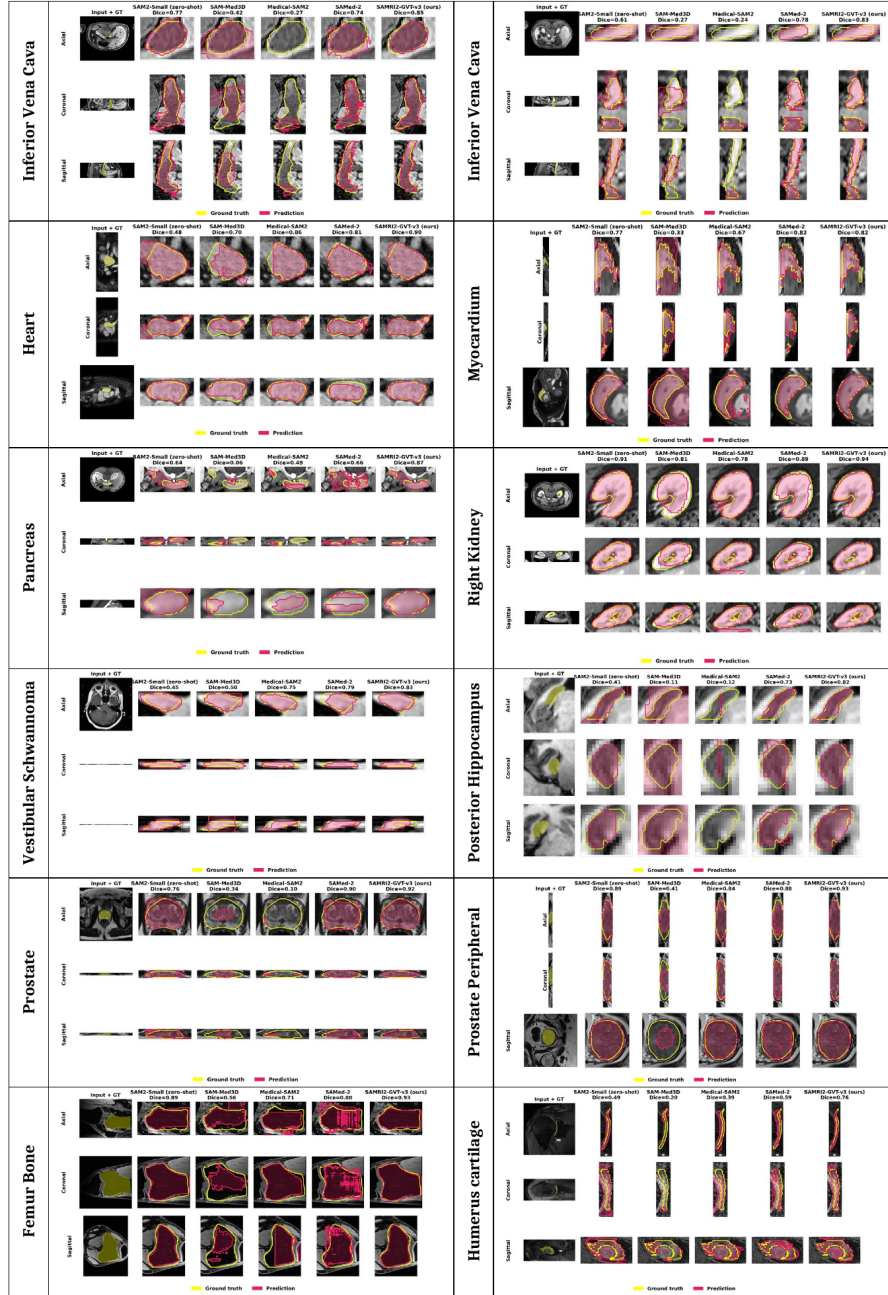


Figure S2. Additional qualitative segmentation examples (Positive).

Figure S2. Additional qualitative segmentation examples (successful cases).

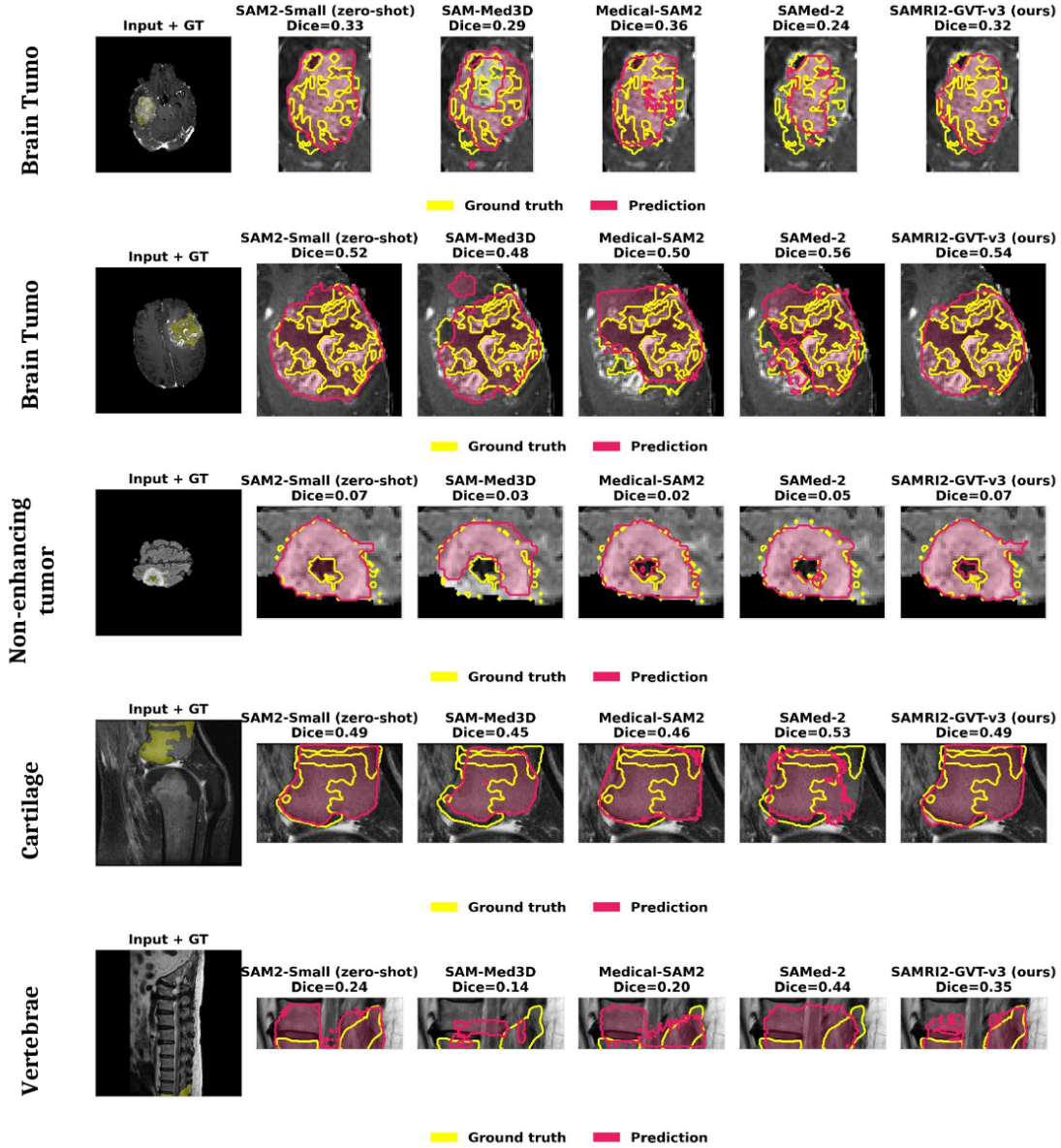


Figure S3. Additional qualitative segmentation examples (Negative).

Figure S3. Additional qualitative segmentation examples (failure cases).

Dataset References

- [1] A. L. Simpson et al., "A large annotated medical image dataset for the development and evaluation of segmentation algorithms," arXiv:1902.09063, 2019; M. Antonelli et al., "The Medical Segmentation Decathlon," Nat. Commun., 13:4128, 2022.
- [2] B. H. Menze et al., "The multimodal brain tumor image segmentation benchmark (BRATS)," IEEE Trans. Med. Imag., 34(10):1993–2024, 2015; S. Bakas et al., "Advancing the cancer genome atlas glioma MRI collections with expert segmentation labels and radiomic features," Sci. Data, 4:170117, 2017; S. Bakas et al., "Segmentation labels and radiomic features for the pre-operative scans of the TCGA-GBM and TCGA-LGG collections," The Cancer Imaging Archive, 2017; S. Bakas et al., "Identifying the best machine learning algorithms for brain tumor segmentation. . .," arXiv:1811.02629, 2018.
- [3] M. R. Hernandez Petzsche et al., "ISLES 2022: A multi-center MRI stroke lesion segmentation dataset," Sci. Data, 9:762, 2022.
- [4] H. J. Kuijff et al., "Standardized assessment of automatic segmentation of white matter hyperintensities and results of the WMH segmentation challenge," IEEE Trans. Med. Imag., 38(11):2556–2568, 2019.
- [5] A. Saha, M. Hosseinzadeh, H. Huisman, "End-to-end prostate cancer detection in bpMRI via 3D CNNs. . .," Med. Image Anal., 73:102155, 2021 (PI-CAI).
- [6] A. S. Becker et al., "Variability of manual segmentation of the prostate in axial T2-weighted MRI: a multi-reader study," Eur. J. Radiol., 121:108716, 2019.
- [7] Prostate (ADC) data from the CVPR 2024 challenge "Segment Anything in Medical Images on Laptop" (MedSAM on Laptop): J. Ma, Y. Zhou, B.

- Wang et al., “Segment Anything in Medical Images and Videos: Benchmark and Deployment,” arXiv:2408.03322, 2024; proceedings in Medical Image Segmentation Foundation Models (CVPR 2024 Challenge), LNCS 15458, Springer, 2025.
- [8] A. Fedorov et al., “An annotated test-retest collection of prostate multiparametric MRI,” *Sci. Data*, 5:180281, 2018; hosted on The Cancer Imaging Archive: K. Clark et al., “The Cancer Imaging Archive (TCIA): maintaining and operating a public information repository,” *J. Digit. Imaging*, 26(6):1045–1057, 2013.
- [9] G. Litjens et al., “Evaluation of prostate segmentation algorithms for MRI: the PROMISE12 challenge,” *Med. Image Anal.*, 18(2):359–373, 2014.
- [10] F. Ambellan et al., “Automated segmentation of knee bone and cartilage combining statistical shape knowledge and CNNs: Data from the Osteoarthritis Initiative,” *Med. Image Anal.*, 52:109–118, 2019; source imaging from the Osteoarthritis Initiative (OAI): C. G. Peterfy, E. Schneider, M. Nevitt, “The osteoarthritis initiative: report on the design rationale for the MRI protocol for the knee,” *Osteoarthritis Cartilage*, 16(12):1433–1441, 2008.
- [11] Osteoarthritis Initiative (OAI): C. G. Peterfy, E. Schneider, M. Nevitt, “The osteoarthritis initiative: report on the design rationale for the MRI protocol for the knee,” *Osteoarthritis Cartilage*, 16(12):1433–1441, 2008; <https://nda.nih.gov/oai>.
- [12] O. Bernard et al., “Deep learning techniques for automatic MRI cardiac multi-structures segmentation and diagnosis: is the problem solved?,” *IEEE Trans. Med. Imag.*, 37(11):2514–2525, 2018.
- [13] D. F. Pace et al., “HVS MR-2.0: A 3D cardiovascular MR dataset for whole-heart segmentation in congenital heart disease,” *Sci. Data*, 11:721, 2024.
- [14] T. A. D’Antonoli et al., “TotalSegmentator MRI: robust sequence-independent segmentation of multiple anatomic structures in MRI,” *Radiology*, 314(2):e241613, 2025.
- [15] R. Dorent et al., “CrossMoDA 2021 challenge: benchmark of cross-modality domain adaptation for vestibular schwannoma and cochlea segmentation,” *Med. Image Anal.*, 83:102628, 2023; original dataset: J. Shapey et al., “Segmentation of vestibular schwannoma from MRI, an open annotated dataset and baseline algorithm,” *Sci. Data*, 8:286, 2021.
- [16] S. S. Chandra et al., “Focused shape models for hip joint segmentation in 3D magnetic resonance images,” *Med. Image Anal.*, 18(3):567–578, 2014.
- [17] D. Zukić et al., “Robust detection and segmentation for diagnosis of vertebral diseases using routine MR images,” *Comput. Graph. Forum*, 33(6):190–204, 2014.
- [18] Y. Ji et al., “AMOS: a large-scale abdominal multi-organ benchmark for versatile medical image segmentation,” *Adv. Neural Inf. Process. Syst.*, 35:36722–36732, 2022.
- [19] S. R. Bowen et al., “Tumor radiomic heterogeneity: multiparametric functional imaging to characterize variability and predict response following cervical cancer radiation therapy,” *J. Magn. Reson. Imaging*, 47(5):1388–1396, 2018.