

Liver Metastasis Detection at Reduced Radiation Dose: Diagnostic Evaluation of a Novel Organ-Level Tube Current Modulation Method

Maria Jose Medrano^{1*}; Sen Wang^{1*}; Kaylee W. Fang¹; Jennie Cao¹; Liyan Sun¹; Luyao Shen¹; Andy Shon¹; Luke Yoon¹; Grant M. Stevens²; Adam S. Wang¹; Justin R. Tse¹

¹Department of Radiology, Stanford University School of Medicine, Stanford University, Stanford, California, USA

²GE HealthCare, Menlo Park, California, USA

*Maria Jose Medrano and Sen Wang contributed equally to this work.

Abstract

Objective: To evaluate a novel organ-level tube current modulation (TCM) method (retroOpt) designed to minimize effective radiation dose while preserving the diagnostic image quality required for accurate liver metastasis detection.

Materials and Methods: In this single-center, institutional review board-approved, retrospective study, 22 patients with 68 liver lesions (38 malignant, 30 benign) underwent portal venous phase contrast-enhanced chest-abdomen-pelvis CT. Five image series were generated per patient using projection-domain noise emulation: 1) original full-dose acquisition (Orig-100), 2) uniform dose reduction at 40% (UD-40) and 3) 60% (UD-60) of the original effective dose, and 4) optimized organ-level TCM at 40% (Opt-40) and 5) 60% (Opt-60) of the original effective dose. Three abdominal radiologists independently detected and classified lesions and assessed image quality using a 5-point Likert scale. Per-lesion sensitivity was compared using the McNemar test. Malignant lesion-size thresholds corresponding to 50% and 90% sensitivity (x_{50} and x_{90}) were estimated using binomial logistic regression. Image quality scores were compared using the Wilcoxon signed-rank test.

Results: At 60% of the original effective dose, organ-level TCM (Opt-60) achieved mean malignant lesion detection sensitivity comparable to full-dose CT (82% vs 80%) and improved malignant lesion detection sensitivity over equivalent uniform dose reduction (67%). At the 40% dose level, organ-level TCM (Opt-40) outperformed uniform dose reduction (UD-40), increasing mean malignant lesion detection sensitivity from 52% to 68%. Similar gains were observed for mean all-lesion detection, improving from 54% (UD-40) to 64% (Opt-40) at the 40% dose

level and from 65% (UD-60) to 77% (Opt-60) at the 60% dose level. Sensitivity for all lesions with Opt-60 (77%) remained comparable to full-dose CT (76%). Optimized-versus-uniform comparisons were statistically significant for malignant and all-lesion detection (all $p \leq 0.003$). Malignant lesion-size thresholds for Opt-60 ($x_{50} = 3.3$ mm; $x_{90} = 15.7$ mm) closely matched those of Orig-100 ($x_{50} = 3.7$ mm; $x_{90} = 16.8$ mm). Median image quality ratings were consistently higher with organ-level TCM than with equivalent uniform dose reduction, and median Opt-60 ratings ranged from 3 to 4.5 across readers.

Conclusions: Organ-level TCM preserved liver metastasis detection and diagnostic image quality at 60% of the original effective dose, outperforming conventional uniform dose reduction. Task-specific organ-level dose optimization may enable greater CT dose reduction than conventional uniform strategies for patients requiring repeated metastases surveillance.

Keywords: Liver metastases, computed tomography, tube current modulation, effective dose, lesion detectability

1. Introduction

The liver represents one of the most common sites of metastases across all cancers.¹⁻³ Contrast-enhanced CT is the primary imaging modality used for their detection and surveillance.⁴ Most hepatic metastases enhance differently relative to surrounding parenchyma, but often the attenuation differences are subtle and make them challenging to detect on CT.^{5,6} Consequently, accurate liver metastasis detection requires sufficient image quality and contrast-to-noise ratio.

Noise can be reduced in CT by increasing radiation dose. However, this represents a trade-off in long-term safety, particularly in oncology patients who are expected to undergo lifelong surveillance CT.⁷ Thus, an ongoing goal has been to reduce radiation exposure while preserving diagnostic image quality in liver metastasis detection. Prior strategies have used conventional filtered back projection, advanced iterative reconstruction, and denoising techniques.⁸⁻¹¹ However, these approaches primarily quantified dose using scanner-output metrics such as CT dose index volume (CTDI_{vol}) and dose-length product (DLP), or effective dose estimated from DLP using population-average conversion factors.¹² Such metrics reflect overall scanner output or generic, population-level risk estimates rather than the radiation burden specific to an individual patient's anatomy. Furthermore, few studies have examined organ-level tube current modulation (TCM) optimized for hepatic lesion detection as a strategy to further reduce patient radiation dose.^{13,14}

A recently proposed organ-level TCM method optimizes tube current profiles to minimize effective dose, an individualized radiation-risk metric computed from organ dose estimates and maximize image quality for selected organs of interest.¹⁵ The objective of this study is to evaluate whether this novel organ-level TCM approach can reduce effective radiation dose but also preserve diagnostic image quality required to accurately detect liver metastases.

2. Materials and Methods

2.1 Patient Selection

This single-center, retrospective study was HIPAA compliant and approved by our institutional review board. From March 13, 2023 to June 23, 2023, 287 patients at a comprehensive cancer center underwent outpatient CT chest-abdomen-pelvis in single-energy mode. A research coordinator then reviewed the electronic medical record

(EMR) and excluded patients who did not have a history of cancer (n=11). The research coordinator reviewed the CT reports of patients with a history of cancer to determine if liver metastases were present (n=38). A board certified, fellowship-trained, abdominal radiologist with 4 years of post-training experience subsequently reviewed the original clinical images and the EMR on a clinical workstation to exclude patients based on the following criteria: 1) widespread liver metastases (defined as >7 metastases) and 2) lack of a reference standard, defined as tissue confirmation or subsequent imaging follow-up (MRI, CT, or PET/CT). For patients who underwent multiple CT exams during the time period, the first CT was selected. At this stage, all liver lesions (benign and malignant) were annotated, and the long-axis dimension of each lesion was measured in the axial plane. After applying all exclusion criteria, this resulted in 11 patients with 38 metastases. A subset of control cases with cancer but without liver metastases was selected to match the age, sex, and body-mass index of the final patient cohort. Figure 1 depicts the process to derive the patient cohort.

2.2 CT Acquisition and Image Reconstruction

Patients were imaged on a GE HealthCare Revolution Apex (Waukesha, WI) scanner using a single-energy CT chest-abdomen-pelvis protocol at our institution. Weight-based iodinated contrast (Isovue-370, Bracco Diagnostics Inc., Monroe Township, NJ, USA) was administered intravenously and images were obtained in the portal venous phase. CT exams were acquired in helical mode using the vendor's automatic exposure control (SmartmA, GE HealthCare), which modulates scanner output based on the scout scan and a weight-based noise index defining the target image quality. Table I summarizes acquisition parameters, including CTDI_{vol}, effective dose, and noise index. Hereafter, the original effective dose corresponding to the scanner's proprietary automatic exposure control is referred to as *Orig-100*. Original reconstructed images, along with projection and tube current modulation data, of the final patient cohort and controlled cases were exported from the clinical CT scanner at our institution.

A previously proposed organ-level tube current modulation optimization method¹⁵, hereafter termed *retroOpt AEC*, was retrospectively implemented to generate optimal modulation profiles for 40% and 60% effective dose (*Opt-40* and *Opt-60*). In this framework, organ-specific weights can be assigned to preferentially allocate dose and thus enhance image quality to organs of interest while maintaining a desired overall effective dose reduction (Supplementary Materials, Figure S1). In this study, organ weights in the effective dose calculation were selected

to ensure that the liver dose matched the minimum levels reported in the literature required to maintain diagnostic image quality for hepatic lesion detection.^{6,10,16,17} Low-dose images corresponding to Opt-40 and Opt-60 were subsequently reconstructed using a projection-domain noise emulation method at the respective reduced dose levels, followed by image reconstruction.¹⁸ All reconstructions were performed by filtered backprojection with a Standard kernel. To evaluate the added benefit of organ-level dose optimization over uniform dose reduction, additional images at 40% and 60% of the original effective dose were reconstructed using the same low-dose emulation method, with dose distribution following the vendor's original TCM (denoted as *UD-40* and *UD-60*). This overall pipeline resulted in five reconstructions per patient for comparative analysis: *Orig-100*, *Opt-40*, *Opt-60*, *UD-40*, and *UD-60*.

2.3 Image Evaluation

Three board-certified, fellowship-trained abdominal radiologists with 7, 8, and 17 years of post-training experience independently reviewed CT images across different dose levels and TCM schemes in 5 sessions. A randomized, blinded reading schedule was generated for each reader to evaluate one exam corresponding to UD-40, UD-60, Opt-40, Opt-60, or Orig-100 per session for each patient. Each session was separated by a minimum 30-day wash-out period to reduce recall bias. The imaging consisted of axial images at 2.5 mm section thickness and reformations in the coronal and sagittal plane at 2 mm section thickness, respectively. All imaging was evaluated on a research instance of our institutional PACS. Images were anonymized, and radiologists were blinded to dose level and patient medical records, including histopathology reports, comparison examinations, and if the patients had a history of cancer.

Radiologists were instructed to annotate all liver lesions in the axial plane, classify them as malignant or benign, and measure size. There was no minimum or maximum allowable lesion size. Readers were also asked to separately assign an image quality grade in the abdomen region using a pre-assigned Likert score, with overall diagnostic image quality graded in a 5-point scale (1 = non-diagnostic image quality, 2 = sub-optimal diagnosis with compromised diagnosis confidence, 3 = acceptable diagnostic image quality, 4 = good diagnostic image quality, 5 = excellent diagnostic image quality). Supplementary material S2 provides further details on the instructions followed by each radiologist to define image quality based on the defined Likert scale.

2.4. Statistical Analysis

The size and location of each lesion detected by the readers were recorded for each patient across all dose levels and TCM schemes. Per-reader sensitivity for malignant lesions and for all lesions (malignant and benign combined) was computed on a per-lesion basis by assigning a binary detection outcome (1 = detected, 0 = not detected) relative to a reference standard established from prior imaging and patient medical records. A liver lesion was considered detected if the reader's lesion annotations overlapped with the corresponding reference. Wilson 95% confidence intervals were computed for all sensitivity estimates. Statistical differences in sensitivity between Opt-40 and UD-40, and between Opt-60 and UD-60, were assessed using McNemar's paired test, with $p < 0.05$ considered statistically significant.

Aggregated sensitivity for malignant lesions and for all lesions was additionally evaluated as a function of lesion size and modeled using a binomial logistic regression. Bootstrap-based 95% confidence intervals were obtained by resampling lesions with replacement, refitting the model for each iteration, and extracting the 2.5th and 97.5th percentiles of the resulting distributions. From the fitted models, lesion size thresholds corresponding to 50% and 90% sensitivity were derived for each condition, along with their bootstrap confidence intervals.

Median Likert scores for abdominal image quality were calculated for each reader across the 22 patient cases. Differences in subjective image quality between UD-40 and Opt-40, and between UD-60 and Opt-60, were evaluated using the Wilcoxon signed-rank test, with $p < 0.05$ indicating statistical significance. All statistical analyses were performed in Python 3.10 using the SciPy 1.15 and Statsmodels 0.14 libraries.

3. Results

3.1 Patient Cohort

The cohort consisted of 68 liver lesions (30 benign, 38 malignant) among 22 patients. The patients had the following primary sites of cancer: colorectal (n=7), lung (n=4), ovarian (n=2), stomach (n=3), breast (n=2), cervical (n=1), endometrial (n=1), lymphoma (n=1), and prostate (n=1). Benign liver lesions included simple cysts, hemangiomas, and calcified granulomas. Overall, patients had a median of 3 liver lesions (range: 0–9) with a median lesion size of 1.8 cm (range: 0.2–9.4 cm). By size, 35 lesions measured <1 cm, 18 lesions measured between 1-2

cm, and 15 lesions measured >2 cm. The patient cohort characteristics and their lesions are summarized in Table II.

3.2 Diagnostic Performance

Detection rates for the readers are summarized in Table III. Malignant lesions were detected at higher percentage rates with Opt-40 (mean = 68, range = 63–74) than with UD-40 (mean = 52, range = 47–55). Similarly, Opt-60 resulted in higher malignant lesion detection rates (mean = 82, range = 68–90) than UD-60 (mean = 67, range = 58–79), and remained comparable to Orig-100 (mean = 80, range = 71–90). A similar trend was observed for the detection of all lesions (malignant and benign combined). Specifically, mean detection rates improved from 54 (range = 47–60) with UD-40 to 64 (range = 59–69) with Opt-40, and from 65 (range = 57–72) with UD-60 to 77 (range = 69–84) with Opt-60. Furthermore, all-lesion detection rates for Opt-60 were comparable to those obtained with full-dose Orig-100 images (mean = 76, range = 72–82). McNemar’s paired test showed statistically significant differences in sensitivity between UD-40 and Opt-40 and between UD-60 and Opt-60; these differences were significant for malignant lesions ($p = 0.0005$ and $p = 0.003$, respectively) and all lesions (i.e. both benign and malignant) ($p = 0.002$ and $p = 0.0002$, respectively).

3.3 Relationship of Diagnostic Performance with Lesion Size and Anatomical Location

Figure 2 provides a visual depiction of the detection rates for the 68 liver lesions (38 malignant and 30 benign) along with their corresponding size, location, and benignity. Sensitivity values for the 8 liver segments and for the different acquisition methods are reported in Supplementary Material S3. Detection rates for the different liver segments improved from UD-40 and Opt-40 to UD-60, Opt-60, and Orig-100, with two exceptions: Opt-60 performed slightly better than Orig-100 in segments II/III and IV, and Opt-40 performed comparably to Opt-60 in segment VII. Otherwise, no appreciable difference was observed among the sensitivity values across the different liver segments.

Figure 3 shows lesion size versus sensitivity plots for malignant lesion detection across the five acquisition methods. Lesion-size thresholds corresponding to 50% and 90% sensitivity (x_{50} and x_{90}) are summarized in Table IV. Smaller malignant lesions were detectable with higher sensitivity in Opt-40 ($x_{50} = 8.7$ mm, $x_{90} = 19.2$ mm) than

in its uniform dose-reduction counterpart, UD-40 (x_{50} = 13.6 mm, x_{90} = 31.7 mm). A similar trend was observed for Opt-60 (x_{50} = 3.3 mm, x_{90} = 15.7 mm) compared with UD-60 (x_{50} = 7.9 mm, x_{90} = 24.0 mm). Furthermore, thresholds were comparable between Opt-60 and Orig-100 (x_{50} = 3.7 mm, x_{90} = 16.8 mm).

A similar trend was observed for the combined sensitivity of all liver lesions (i.e., both benign and malignant lesions). As shown on Table IV, smaller lesions were detectable with higher sensitivity in Opt-40 (x_{50} = 8.1 mm, x_{90} = 17.8 mm) than in its uniform dose reduction counterpart, UD-40 (x_{50} = 11.1 mm, x_{90} = 25.4 mm). This can also be observed for Opt-60 (x_{50} = 4.2 mm, x_{90} = 13.7 mm) compared with UD-60 (x_{50} = 7.5 mm, x_{90} = 18.4 mm). Similarly, lesion size thresholds corresponding to 50% and 90% sensitivity were comparable between Opt-60 (x_{50} = 4.2 mm, x_{90} = 13.7 mm) and Orig-100 (x_{50} = 3.6 mm, x_{90} = 15.7 mm). Supplementary Material S4 shows lesion size versus sensitivity plots for detection of all lesions (benign and malignant) across all five acquisition methods.

3.4. Image Quality

The median image quality ratings for UD-40 across the three readers ranged from suboptimal with compromised diagnostic confidence, (Likert Score = 2) to acceptable diagnostic quality (Likert Score = 3). In contrast, median ratings for Opt-40 ranged from diagnostic with acceptable image quality (Likert Score = 3) to good image quality (Likert Score = 4). UD-60 images were consistently rated by all readers as diagnostic with acceptable image quality (Likert score = 3). Opt-60 images received higher quality ratings, ranging from acceptable to good diagnostic image quality (Likert score = 3–4.5). Orig-100 received the best median image quality ratings, ranging from good to excellent (Likert Score = 4–5). Figure 4 summarizes the image quality ratings assigned by the readers across the different acquisition methods. Based on the Wilcoxon Signed-Rank Test, differences in median Likert Score ratings between UD-40 and Opt-40 were statistically significant for Readers 1 and 3 ($p \leq 0.012$) but not for Reader 2 ($p = 0.24$). Differences between UD-60 and Opt-60 were statistically significant for all readers ($p \leq 0.0067$).

Figures 5 and 6 present representative cases illustrating the improvement in image quality from uniform (UD) to optimized (Opt) dose reduction across the UD-40, Opt-40, UD-60, Opt-60, and Orig-100 series. In both figures, image quality improved in that order, as did lesion detection.

4. Discussion

We demonstrate that retroOpt, a retrospective organ-level tube current modulation algorithm, can maintain hepatic lesion detectability at 60% of the original effective dose, or a 40% reduction in effective dose. Unlike conventional TCM methods that primarily adjust tube current according to patient attenuation, retroOpt performs task-specific optimization to preserve image quality within organs of interest while maintaining a predefined effective dose target.¹⁵ The lesion detectability among our 3 readers was comparable to full-dose acquisition and outperformed equivalent uniform dose reduction strategies. This pattern was also observed in our subjective image quality results, where optimized acquisitions were rated higher than their uniform dose-reduction counterparts. The benefit was most pronounced for smaller lesions, with optimization reducing the lesion-size thresholds required for 50% and 90% sensitivity relative to uniform dose reduction. For Opt-60, image quality and detection thresholds also approached those observed at full dose.

This ability to reduce dose while preserving detectability is especially important in liver imaging. The liver is among the most common organs of metastatic spread, and metastatic disease accounts for about 90% of cancer-related deaths.^{1-3,19} Consequently, repeated contrast-enhanced CTs are routinely performed for detection and surveillance of hepatic metastases. However, because liver lesions frequently demonstrate only subtle attenuation differences relative to surrounding parenchyma, detection is particularly vulnerable to increased image noise at reduced dose levels.

Prior studies have attempted to define dose boundaries for which reliable detection remains feasible. Filtered back projection reconstruction techniques demonstrated that liver metastasis detection is largely preserved at moderate dose reductions of 25–50% (CTDI_{vol}: 8.1–10.4 mGy) but deteriorates substantially at more aggressive reductions corresponding to CTDI_{vol} levels of 1.7–5.0 mGy.^{6,16,17,20-22} Our results with uniform dose reduction corroborate this trend. Relative to Orig-100, sensitivity loss at UD-60 (CTDI_{vol} = 7.7 mGy) remained modest, whereas substantially greater degradation was observed at UD-40 (CTDI_{vol} = 5.2 mGy). This behavior is consistent with findings from Fält et al., who reported progressively larger sensitivity losses (10%-20%) below a CTDI_{vol} threshold of approximately 8.6 mGy.²³

RetroOpt may effectively shift this threshold relative to conventional uniform dose reduction strategies which would enable preservation of lesion detectability at lower dose levels. Detection sensitivity for Opt-40 differed by

only 1% from that of UD-60 despite acquisition at a 20% lower overall effective dose. Similarly, sensitivity differences between Opt-60 and Orig-100 were minimal (2%); lesion detectability was largely preserved despite a 40% dose reduction. The magnitude of this difference was comparable to that previously reported between 75% and full dose in prior studies^{6,16,23}. Clinically, these results suggest that retroOpt effectively recovered an additional 15% dose reduction while incurring a similar loss in detectability.

Smaller lesions were more sensitive to image quality degradation. Fletcher et al. previously reported reduced reader confidence for malignant hepatic lesions below 15 mm.^{10,16} Nevertheless, in our study, the diameter thresholds corresponding to 50% and 90% sensitivity were comparable between Opt-60 ($x_{50} = 3.3$ mm, $x_{90} = 15.7$ mm) and Orig-100 ($x_{50} = 3.7$ mm, $x_{90} = 16.8$ mm). Logistic fits further revealed improved sensitivity for lesions in the 5–15 mm range with Opt-60 relative to all other dose-reduction conditions. Our results indicate that our organ-level optimization can aid the detection of smaller lesions. Consistent with prior work²³, we found no apparent correlation between lesion detectability and location across the different dose reduction strategies. This suggests that improvements in contrast-to-noise ratio and dose reduction strategies were spatially uniform across hepatic segments.

Previous studies of organ-level TCM demonstrated that noise levels and subjective image quality within the liver and adjacent organs could be preserved at moderately reduced overall dose^{13,14}. However, these studies did not evaluate whether these improvements translate into equivalent clinical task performance. The present findings address this gap: retroOpt achieved acceptable-to-good diagnostic image quality in the liver at Opt-60 (Likert scale 3–4.5), and these image quality improvements translated into preserved lesion detection. By preferentially increasing dose deposition within the liver, retroOpt improved the contrast-to-noise conditions for lesion-background discrimination, which represents the primary determinant of detectability for low-contrast hepatic lesions.^{24–26}

Besides evaluating diagnostic task performance, our study also differs from prior work by addressing dosimetric optimization. While previous organ-level TCM approaches optimized dose with respect to $CTDI_{vol}$, retroOpt minimizes effective dose, a unified radiation-risk metric that accounts for the differential biological sensitivity of irradiated tissues.¹² Although retroOpt increased local dose deposition in the liver to preserve diagnostic image

quality, this redistribution produced a net reduction in effective dose because the liver carries a relatively low tissue-weighting factor (0.04) compared to neighboring organs such as the stomach (0.12) and colon (0.12), which are generally not sites of distant organ metastases.²⁷

Our study had several limitations. First, it was a single-center, retrospective study performed on a single CT scanner platform. Generalizability will require validation in larger, multi-center cohorts. Second, the cohort was relatively small (22 patients, 68 lesions). Because each patient required multiple reconstruction conditions, analytical complexity constrained enrollment. The primary goals, however, were to determine whether organ-level dose optimization through retroOpt could preserve lesion detectability relative to uniform dose reduction and to identify the dose threshold at which this benefit was most pronounced, rather than to establish definitive clinical performance estimates. Third, dose reduction was simulated using projection-domain noise insertion rather than prospectively acquired low-dose images to enable comparison across dose levels without additional patient exposure. Our noise model and the retroOpt optimization have been previously validated against phantom acquisitions^{15,28}, preliminarily supporting the accuracy of our results. Fourth, evaluation was limited to filtered back projection. Advanced iterative and deep learning–based reconstructions have previously shown promise at reduced dose^{10,11}, with benefits greatest near the dose threshold at which observer performance begins to degrade.¹⁶ The thresholds identified here will provide a principled basis for evaluating retroOpt in combination with such methods. Finally, the current retroOpt implementation is retrospective. A prospective organ-level modulation would require real-time organ localization and dose estimation. Recent deep learning advances suggest this is feasible^{29–31}; a prospective design would be an important next step toward clinical translation.

In summary, the proposed organ-level TCM approach reduced effective radiation dose while preserving the diagnostic image quality required for accurate detection of hepatic metastases. Furthermore, organ-level TCM optimization may extend dose reduction beyond the threshold achievable with conventional uniform dose reduction.

5. Acknowledgments

This work was supported by GE HealthCare and Stanford Propel Postdoctoral Fellowship.

6. References

1. Hess KR, Varadhachary GR, Taylor SH, et al. Metastatic patterns in adenocarcinoma. *Cancer*. 2006;106(7):1624-1633.
2. Budczies J, von Winterfeld M, Klauschen F, et al. The landscape of metastatic progression patterns across major human cancers. *Oncotarget*. 2015;6(1):570-583.
3. Mielgo A, Schmid MC. Liver Tropism in Cancer: The Hepatic Metastatic Niche. *Cold Spring Harb Perspect Med*. 2020;10(3):a037259.
4. Kaur H, Hindman NM, Al-Refaie WB, et al. ACR Appropriateness Criteria® Suspected Liver Metastases. *Journal of the American College of Radiology*. 2017;14(5):S314-S325.
5. Robinson PJ. Imaging liver metastases: current limitations and future prospects. *Br J Radiol*. 2000;73(867):234-241.
6. Kanal KM, Chung JH, Wang J, et al. Impact of Incremental Increase in CT Image Noise on Detection of Low-contrast Hypodense Liver Lesions. *Acad Radiol*. 2014;21(10):1233-1239.
7. Smith-Bindman R, Chu PW, Azman Firdaus H, et al. Projected Lifetime Cancer Risks From Current Computed Tomography Imaging. *JAMA Intern Med*. 2025;185(6):710.
8. Prakash P, Kalra MK, Kambadakone AK, et al. Reducing Abdominal CT Radiation Dose With Adaptive Statistical Iterative Reconstruction Technique. *Invest Radiol*. 2010;45(4):202-210.
9. Sagara Y, Hara AK, Pavlicek W, Silva AC, Paden RG, Wu Q. Abdominal CT: Comparison of Low-Dose CT With Adaptive Statistical Iterative Reconstruction and Routine-Dose CT With Filtered Back Projection in 53 Patients. *American Journal of Roentgenology*. 2010;195(3):713-719.
10. Fletcher JG, Yu L, Li Z, et al. Observer performance in the detection and classification of malignant hepatic nodules and masses with CT image-space denoising and iterative reconstruction. *Radiology*. 2015;276(2):465-478.
11. Kanan A, Pereira B, Hordonneau C, et al. Deep learning CT reconstruction improves liver metastases detection. *Insights Imaging*. 2024;15(1).
12. McCollough C, Cody D, Edyvean S, et al. *AAPM Report No. 96: The Measurement, Reporting, and Management of Radiation Dose in CT*. 2008.
13. Devapalasundaram A, Lau KK, Paul E, et al. Liver detection algorithm: its efficacy for CT noise reduction in the liver. *Abdominal Radiology*. 2016;41(3):493-499.
14. Fillon M, Si-Mohamed S, Coulon P, Vuillod A, Klahr P, Boussel L. Reduction of patient radiation dose with a new organ based dose modulation technique for thoraco-abdominopelvic computed tomography (CT) (Liver dose right index). *Diagn Interv Imaging*. 2018;99(7-8):483-492.
15. Wang S, Medrano Matamoros MJ, Imran AAZ, Stevens G, Tse JR, Wang AS. Retrospective tube current modulation optimization of individualized organ-level CT dose and image quality. In: Fahrig R, Sabol JM, Li K, eds. *Medical Imaging 2024: Physics of Medical Imaging*. SPIE; 2024:41.
16. Fletcher JG, Yu L, Fidler JL, et al. Estimation of Observer Performance for Reduced Radiation Dose Levels in CT. *Acad Radiol*. 2017;24(7):876-890.
17. Wessling J, Esseling R, Raupach R, et al. The effect of dose reduction and feasibility of edge-preserving noise reduction on the detection of liver lesions using MSCT. *Eur Radiol*. 2007;17(7):1885-1891.

18. Benson TM, De Man BKB. Synthetic CT noise emulation in the raw data domain. In: *IEEE Nuclear Science Symposium & Medical Imaging Conference*. IEEE; 2010:3169-3171.
19. Chaffer CL, Weinberg RA. A Perspective on Cancer Cell Metastasis. *Science* (1979). 2011;331(6024):1559-1564.
20. Pooler BD, Lubner MG, Kim DH, et al. Prospective Evaluation of Reduced Dose Computed Tomography for the Detection of Low-Contrast Liver Lesions: Direct Comparison with Concurrent Standard Dose Imaging. *Eur Radiol*. 2017;27(5):2055-2066.
21. Son JH, Kim SH, Yoon JH, Lee Y, Lim YJ, Kim SJ. Comparison of Model-Based Iterative Reconstruction, Adaptive Statistical Iterative Reconstruction, and Filtered Back Projection for Detecting Hepatic Metastases on Submillisievert Low-Dose Computed Tomography. *J Comput Assist Tomogr*. 2017;41(4):644-650.
22. Pickhardt PJ, Lubner MG, Kim DH, et al. Abdominal CT With Model-Based Iterative Reconstruction (MBIR): Initial Results of a Prospective Trial Comparing Ultralow-Dose With Standard-Dose Imaging. *American Journal of Roentgenology*. 2012;199(6):1266-1274.
23. Fält T, Söderberg M, Hörberg L, et al. Simulated dose reduction for abdominal CT with filtered back projection technique: Effect on liver lesion detection and characterization. *American Journal of Roentgenology*. 2019;212(1):84-93.
24. Jang HJ, Kim TK, Wilson SR. Imaging of malignant liver masses: characterization and detection. *Ultrasound Q*. 2006;22(1):19-29.
25. Garrett R. Solid Liver Masses: Approach to Management from the Standpoint of a Radiologist. *Curr Gastroenterol Rep*. 2013;15(12):359.
26. Vialle R, Boucebei S, Richer JP, et al. Preoperative detection of hepatic metastases from colorectal cancer: Prospective comparison of contrast-enhanced ultrasound and multidetector-row computed tomography (MDCT). *Diagn Interv Imaging*. 2016;97(9):851-855.
27. Valentin J. ICRP Publication 103: The 2007 recommendations of the international commission on radiological protection. *Ann ICRP*. 2007;37(2-4):1-133.
28. Wang S, Medrano MJ, Imran AAZ, et al. Automated estimation of individualized organ-specific dose and noise from clinical CT scans. *Phys Med Biol*. 2025;70(3):035014.
29. Medrano MJ, Wang S, Sun L, et al. Scout-Dose-TCM: direct and prospective scout-based estimation of personalized organ and effective doses from tube current modulated CT exams. *Phys Med Biol*. 2025;70(24):245011.
30. Imran AAZ, Wang S, Pal D, et al. Scout-Net: Prospective Personalized Estimation of CT Organ Doses from Scout Views. *arXiv:231215354*. Published online 2023.
31. Sun L, Medrano MJ, Fang KW, Stevens GM, Tse JR, Wang AS. Deep Learning Enhancement of Ultralow-Dose 3D CT Scouts for Organ Localization and Dose. In: *The 9th International Conference on Image Formation in X-Ray Computed Tomography*. 2026.

6. Tables

Table I: Acquisition and reconstruction parameters.

Acquisition Parameters	Values
Tube voltage (number of patients)	
100 kVp	5 (22.7%)
120 kVp	16 (72.7%)
140 kVp	1 (4.5%)
Median Noise Index (IQR)	12.4 (11.5–13.0)
Median CTDI _{vol} (IQR)	12.9 (10.1–16.3) mGy
Median Effective Dose (IQR)	13.1 (11.0–14.1) mSv
Field-of-View	500 mm
Section Thickness	2.5 mm
Matrix Size	512x512 pixels
Reconstruction Kernel	Standard
Reconstruction	Filtered Backprojection

Table II: Patient baseline characteristics. Values report absolute number with interquartile range or percentages in parentheses.

Parameters	Values
<u>Patients</u>	
Number of patients	22
Median Age (IQR)	57 (43–65) years
<u>Sex</u>	
Male	10 (45%)
Female	12 (55%)
Median BMI (IQR)	25.8 (22.0-29.3) kg/m ²
<u>Cancer Types</u>	
-Colorectal	7 (31.8%)
-Lung	4 (18.2 %)
-Stomach	3 (13.6%)
-Breast	2 (9.1%)
-Ovarian	2 (9.1%)
-Cervical	1 (4.5%)
-Endometrial	1 (4.5%)
-Lymphoma	1 (4.5%)
-Prostate	1 (4.5%)
<u>Lesions</u>	
Number of lesions	68
Malignant	38 (55.9%)
Benign	30 (44.1%)
<u>Number of lesions per patient</u>	
0	2 (9.1%)
1-2	7 (31.8 %)
3-5	9 (40.9%)
>5	4 (18.2%)
Median Lesion Size (IQR)	1.0 (0.7–1.8) cm
<u>Lesion size distribution</u>	
<1 cm	31 (45.6%)
1-2 cm	21 (30.9%)
>2 cm	16 (23.5%)

Table III: Sensitivity performance for different dose reduction methods and reduction levels.

Mean Sensitivity for Malignant Lesions [Confidence Intervals] (%)					
	UD-40	Opt-40	UD-60	Opt-60	Orig-100
R1	47 [31 –64]	63 [46 –78]	58 [41 –74]	68 [51 –83]	71 [54 –85]
R2	55 [38 –71]	66 [49 –80]	63 [46 –78]	87 [72 –96]	79 [63 –90]
R3	53 [36 –69]	74 [57 –87]	79 [63 –90]	90 [75 –97]	90 [75 –97]
Mean Sensitivity for All Lesions [Confidence Intervals] (%)					
	UD-40	Opt-40	UD-60	Opt-60	Orig-100
R1	47 [35 –60]	59 [46 –71]	57 [45 –69]	69 [57 –80]	72 [60 –82]
R2	60 [48 –72]	63 [51 –75]	66 [54 –77]	79 [68 –88]	75 [63 –85]
R3	55 [42 –67]	69 [57 –80]	72 [60 –82]	84 [73 –92]	82 [71 –91]

Table IV: Lesion-size thresholds and corresponding confidence intervals for malignant and all lesions.

Sensitivity Size for Malignant Lesions (mm) [Confidence Intervals]					
	UD-40	Opt-40	UD-60	Opt-60	Orig-100
50%	13.6 [7.22, 17.8]	8.65 [4.68, 11.2]	7.91 [-0.86, 11.4]	3.33 [-8.0, 7.12]	3.69 [-15.8, 7.84]
90%	31.7 [23.0, 44.3]	19.2 [13.4, 28.0]	24.0 [15.2, 44.5]	15.7 [10.4, 23.4]	16.8 [9.95, 29.0]
Sensitivity Size for All Lesions (mm) [Confidence Intervals]					
	UD-40	Opt-40	UD-60	Opt-60	Orig-100
50%	11.1 [7.95, 14.0]	8.05 [5.85, 9.83]	7.53 [4.59, 9.54]	4.17 [-0.61, 6.07]	3.62 [-3.78, 6.32]
90%	25.4 [18.3, 32.1]	17.8 [13.0, 24.3]	18.4 [13.2, 26.9]	13.7 [9.90, 20.2]	15.7 [10.4, 23.8]

7. Figures Captions

FIGURE 1: Flowchart of the patient selection process, with inclusion and exclusion criteria.

FIGURE 2: Schematic of the anatomic locations of all malignant and benign liver lesions across all five acquisition methods, with corresponding lesion type and size. Grayscale shading denotes the number of readers who detected each lesion. Lesions were mapped to their corresponding location defined as the caudate lobe (segment I), left lateral sector (segments II-III), left medial sector (segment IV), right anterior sector (segments V and VIII), and right posterior sector (segments VI and VII). Lesions annotated as adjacent to the capsule or fissure were within 1 cm of either structure.


FIGURE 3: Lesion size versus sensitivity plots for malignant lesion detection across all five acquisition methods. Blue dots denote sensitivity per lesion, the blue line corresponds to the binomial logistic regression fit, and the orange overlay represents the 95% confidence interval. Vertical dotted lines indicate the lesion sizes x_{50} and x_{90} , corresponding to malignant lesion sensitivities of 50% and 90%, respectively. Negative lower bounds on the confidence intervals are a mathematical artifact of the logistic regression.

FIGURE 4: Violin plot summarizing the readers' image quality scores for each acquisition method. The white-filled circle denotes the median score, and the width and number of dots represent the number of responses for each image quality score.

FIGURE 5: Representative images of a malignant liver lesion (solid white arrow), emulated with (A) UD-40, (B) Opt-40, (C) UD-60, (D) Opt-60, and (E) Orig-100. The liver lesion was detected by no readers on UD-40 or Opt-40, by 2 readers on UD-60, by all 3 readers on Opt-60, and by 2 readers on Orig-100. Median Likert image quality scores of 3, 4, 3, 4, and 4 were assigned to panels (A)–(E), respectively.

FIGURE 6: Representative images of two differently sized malignant liver lesions (2.8-cm lesion, solid black arrow; 1.2-cm lesion, solid white arrow), emulated with (A) UD-40, (B) Opt-40, (C) UD-60, (D) Opt-60, and (E) Orig-100. A zoomed-in view of the 1.2-cm lesion is shown in the bottom-right inset (yellow outline). The 2.8-cm lesion was detected by 2 readers on UD-40 and all 3 readers on Opt-40, UD-60, Opt-60, and Orig-100 images. The 1.2-cm lesion was detected by no readers on UD-40, by 2 readers on Opt-40, by all three readers on UD-60 and Opt-

60, and by 2 readers on Orig-100. Median Likert image quality scores of 2, 3, 3, 4, and 5 were assigned to panels (A)–(E), respectively.

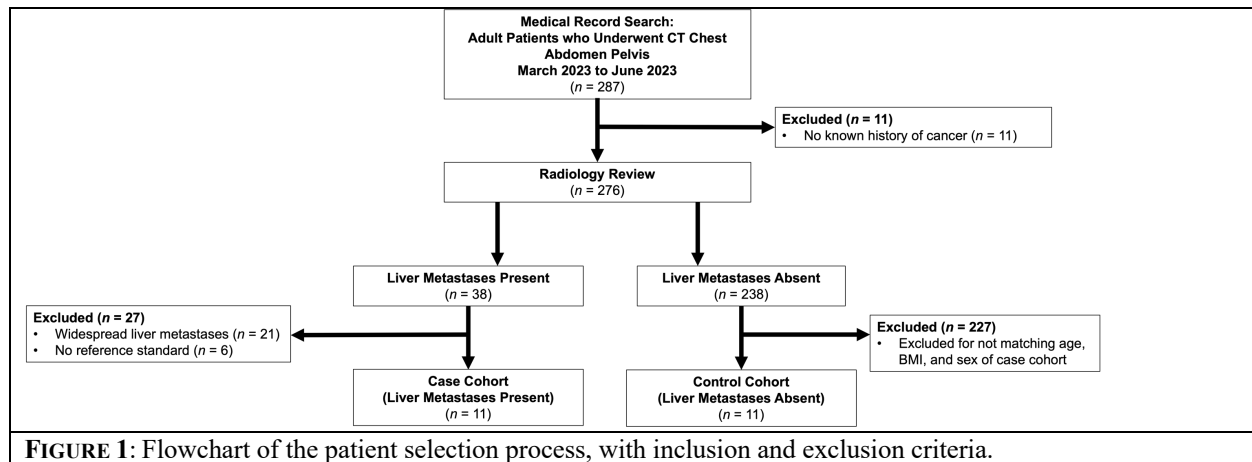

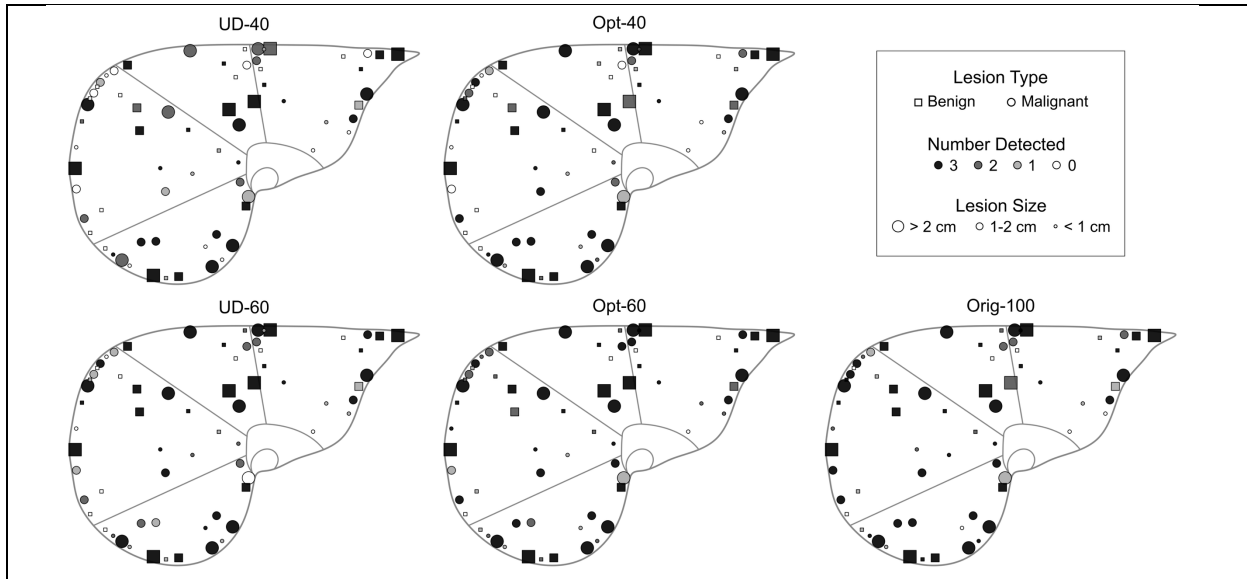


FIGURE 2: Schematic of the anatomic locations of all malignant and benign liver lesions across all five acquisition methods, with corresponding lesion type and size. Grayscale shading denotes the number of readers who detected each lesion. Lesions were mapped to their corresponding location defined as the caudate lobe (segment I), left lateral sector (segments II-III), left medial sector (segment IV), right anterior sector (segments V and VIII), and right posterior sector (segments VI and VII). Lesions annotated as adjacent to the capsule or fissure were within 1 cm of either structure.

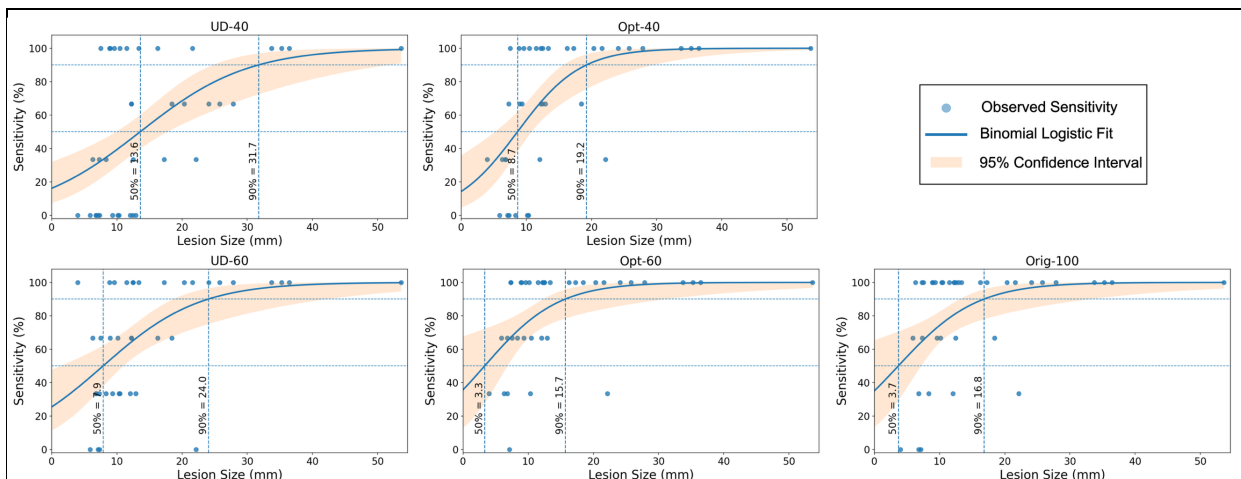


FIGURE 3: Lesion size versus detection sensitivity plots for malignant lesion detection across all five acquisition methods. Blue dots denote sensitivity per lesion, the blue line corresponds to the binomial logistic regression fit, and the orange overlay represents the 95% confidence interval. Vertical dotted lines indicate the lesion sizes x_{50} and x_{90} , corresponding to malignant lesion sensitivities of 50% and 90%, respectively. Negative lower bounds on the confidence intervals are a mathematical artifact of the logistic regression.

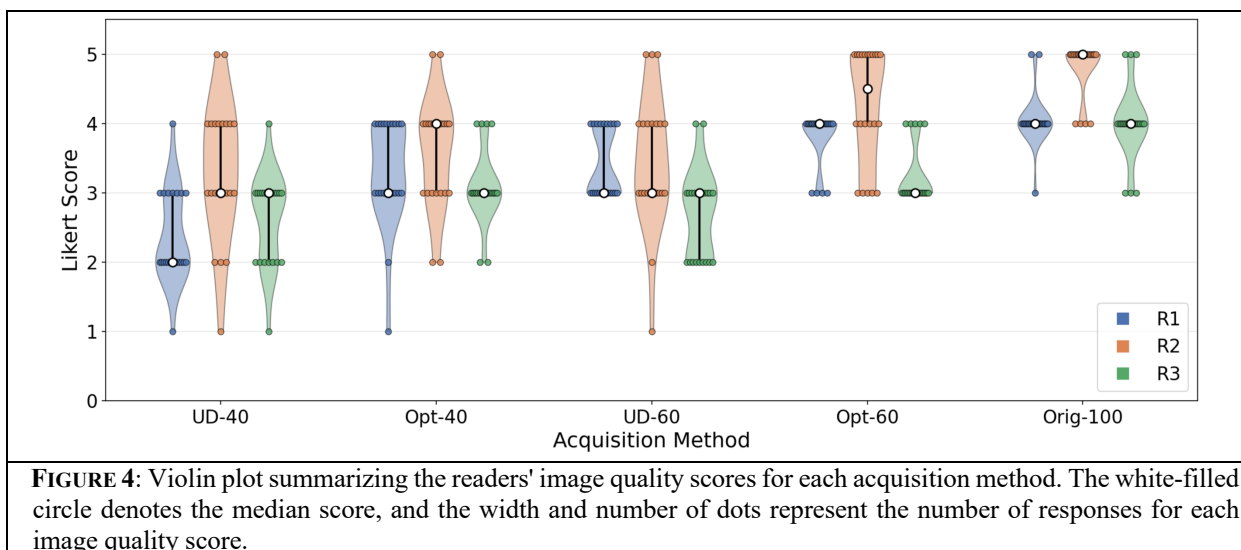


FIGURE 4: Violin plot summarizing the readers' image quality scores for each acquisition method. The white-filled circle denotes the median score, and the width and number of dots represent the number of responses for each image quality score.



FIGURE 5: Representative images of a malignant liver lesion (solid white arrow), emulated with (A) UD-40, (B) Opt-40, (C) UD-60, (D) Opt-60, and (E) Orig-100. The liver lesion was detected by no readers on UD-40 or Opt-40, by 2 readers on UD-60, by all 3 readers on Opt-60, and by 2 readers on Orig-100. Median Likert image quality scores of 3, 4, 3, 4, and 4 were assigned to panels (A)–(E), respectively.


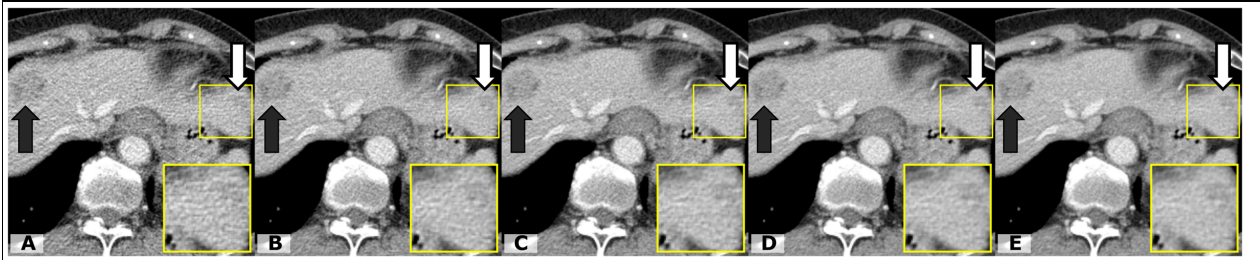


FIGURE 6: Representative images of two differently sized malignant liver lesions (2.8-cm lesion, solid black arrow; 1.2-cm lesion, solid white arrow), emulated with (A) UD-40, (B) Opt-40, (C) UD-60, (D) Opt-60, and (E) Orig-100. A zoomed-in view of the 1.2-cm lesion is shown in the bottom-right inset (yellow outline). The 2.8-cm lesion was detected by 2 readers on UD-40 and all 3 readers on Opt-40, UD-60, Opt-60, and Orig-100 images. The 1.2-cm lesion was detected by no readers on UD-40, by 2 readers on Opt-40, by all three readers on UD-60 and Opt-60, and by 2 readers on Orig-100. Median Likert image quality scores of 2, 3, 3, 4, and 5 were assigned to panels (A)–(E), respectively.

Supplementary Material

S1. Retrospective Organ-Level Tube Current Modulation

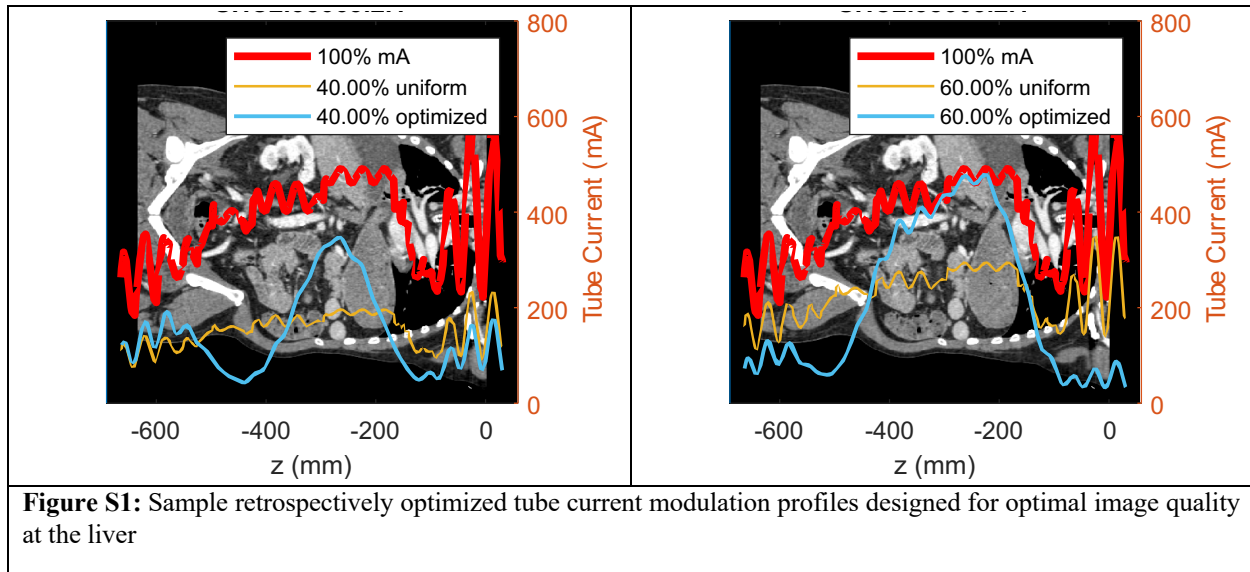

S2. Likert Score Details

Table S1: Detailed description of Likert Score ratings.

5	Diagnostic – Excellent Image Quality - Excellent overall image quality - Very low image noise and reader highly preferred image texture - Very high resolution to depict details with very high confidence. - Very good contrast to differentiate between materials with very high confidence. - Clean images with almost zero image artifacts
4	Diagnostic – Good Image Quality - Good overall image quality - Low image noise and reader preferred image texture - High resolution to depict details with confidence. - Good contrast to differentiate between materials with confidence. - Minor artifacts not interfering with diagnostic decision-making.
3	Diagnostic – Acceptable Image Quality - Acceptable overall image quality - Moderate image noise and acceptable image texture - Acceptable resolution to depict details with acceptable level of confidence. - Acceptable contrast to differentiate between materials with an acceptable level of confidence - Moderate artifacts possibly interfering with diagnostic decision-making.
2	Sub-optimal diagnostic with compromised diagnosis confidence - Poor overall image quality - Severe image noise and poor image texture that can affect diagnosis. - Poor resolution to compromise confidence in depicting details. - Poor contrast to differentiate between materials.

	- Major artifacts affecting visualization of major structures (diagnosis still possible)
1	Non-diagnostic - Unacceptable overall image quality - Unacceptably severe image noise and unacceptable image texture that affects diagnosis. - Unacceptably poor resolution - Unacceptable contrast to differentiate between materials. - Unacceptable artifacts clearly affecting diagnosis. (non-diagnosis)

S3. Sensitivity for Different Lesion Locations

Table S2: Detection sensitivity of benign and malignant lesions across liver segments. No lesions were in the caudate lobe (segment I); detection across segments II–VIII is grouped into anatomic sectors.

Mean Sensitivity for All Lesions Across Different Liver Segments (%)					
Segment	UD-40	Opt-40	UD-60	Opt-60	Orig-100
II/III	56 [41–69]	65 [51–77]	67 [53–79]	78 [64–88]	69 [54–81]
IV a/b	63 [42–81]	59 [39–78]	82 [62–94]	85 [66–96]	78 [58–91]
V/VIII	42 [30–54]	56 [43–67]	57 [45–69]	75 [63–85]	78 [66–87]
VI/VII	65 [50–78]	77 [63–87]	67 [52–79]	77 [63–87]	82 [69–92]

S4. Sensitivity vs. Lesions Size – All Lesions

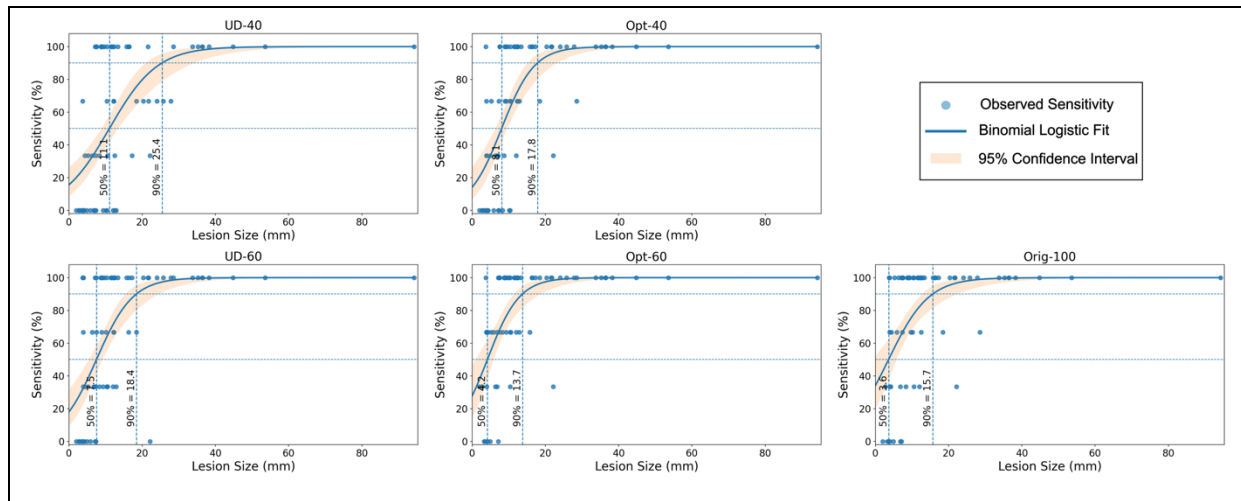


Figure S2: Lesion size versus sensitivity plots for all lesion detection (malignant and benign) across all five acquisition methods. Blue dots denote sensitivity per lesion, the blue line corresponds to the binomial logistic regression fit, and the orange overlay represents the 95% confidence interval. Vertical dotted lines indicate the lesion sizes x_{50} and x_{90} , corresponding to all-lesion sensitivities of 50% and 90%, respectively. Negative lower bounds on the confidence intervals are a mathematical artifact of the logistic regression.